

Instructions for Use of Albumin (ALB) Kit (Bromocresol Green Method)

Package Specification

REF	Reagent	Systems
01.09.00.04.EC.01	R 30 mL x 6	Zybio EXC200/220
01.09.00.04.EC.03	R 60 mL x 2	Hitachi 7180 Zybio EXC400/420

Intended Use

In vitro test for the quantitative determination of albumin (ALB) concentration in human samples (serum). Clinically, it is mainly used as an aid to evaluation of liver function as well as nutritional assessment.

Summary

Albumin is a carbohydrate-free protein, which constitutes 55 - 65% of total plasma protein. It maintains plasma oncotic pressure, and is also involved in the transport and storage of a wide variety of ligands and is a source of endogenous amino acids. Albumin binds and solubilizes various compounds, e.g. bilirubin, calcium and long-chain fatty acids. Furthermore, albumin is capable of binding toxic heavy metal ions as well as numerous pharmaceuticals, which is the reason why lower albumin concentrations in blood have a significant effect on pharmacokinetics.

Hyperalbuminemia is of little diagnostic significance except in the case of dehydration. Hypoalbuminemia occurs during many illnesses and is caused by several factors: compromised synthesis due either to liver disease or as a consequence of reduced protein uptake; elevated catabolism due to tissue damage (severe burns) or inflammation; malabsorption of amino acids (Crohn's disease); proteinuria as a consequence of nephrotic syndrome; protein loss via the stool (neoplastic disease). In severe cases of hypoalbuminemia, the maximum albumin concentration of plasma is 2.5 g/dL (380 µmol/L). Due to the low osmotic pressure of the plasma, water permeates through blood capillaries into tissue (edema). The determination of albumin allows monitoring of a controlled patient dietary supplementation and serves also as an excellent test of liver function.

Principle

Albumin in serum binds to bromocresol green to form a blue-green complex at pH 4.2, which has an absorption peak at the wavelength of 630 nm, and the change in color intensity is directly proportional to the albumin concentration. The albumin concentration in the serum can be obtained by comparing with that in calibrator treated in the same manner.

Reagents Components and Concentration

Components	Main Constituents	Concentration
R	Bromocresol Green	0.15 mmol/L
	Succinic Acid buffer	74.9 mmol/L

The components in different batches are non-interchangeable.

Storage and Validity

- The reagents should be stored at 2 - 8 °C and kept away from direct light and freezing. The unopened reagents are valid for 12 months.
- Once opened, the reagents are stable for 30 days at 2 - 8 °C. For reagents not in use, the cap should be tightened to avoid contamination.
- The production date and expiration date are available on package insert.

System Information

Hitachi 7180, Zybio EXC400/420, Zybio EXC200/220 Chemistry Analyzer. Other models can be used after verification.

Specimen Information

Non-hemolytic serum is suitable for samples, which are stable at 2 - 8 °C for 14 days.

Warnings and Precautions

- For in vitro diagnostic use only. In case of contact, immediately flush the affected area with plenty of water.
- The dosage of reagents and samples can be appropriately increased or decreased according to various instruments, under the condition that the volume proportion of reagents and samples is invariable.
- The necessary precautions should be taken to use the reagents. The bottle cap should be tightened immediately after use to avoid contamination.
- When reagent becomes turbid or the blank absorbance > 0.500, the reagent is failed and should be discarded.
- All samples and reaction wastes should be treated as sources of infection. The waste liquid generated during the experiment and the used packaging materials shall be collected and disposed in accordance with relevant regulations.
- The same sample tested with reagents from different manufacturers may lead to different measured values.
- Warning: This product contains ingredients of human or animal origin. At present, there is no way to completely guarantee that it is free from infectious substances and may be contaminated during use. This product and samples shall be considered as potential sources of infection, and operators should take protective measures and comply with laboratory safety regulations. All waste shall be disposed of in accordance with local regulations.

Test Process

1. Parameters

Method	End-Point Method	Sample/Reagent	1/100
Main Wavelength	630 nm	Reaction Temperature	37 °C
Sub Wavelength	700 nm	Reaction Time	2 min
Reaction Direction	+		

2. Operation

Addition	Blank	Calibration	Detection
Sample (µL)	/	/	3
Calibrator (µL)	/	3	/
Purified Water (µL)	3	/	/
Reagent (µL)	300	300	300
Mix well, measure absorbance A after 2 min.			

3. Calibration

Use Zybio Clinical Chemistry Multi-analyte Calibrator or Randox multi-analyte calibrator. Calibration cycle: Re-calibration is needed in case of replacement of the reagent batch, deviation of quality control, or the major repair and maintenance of equipment.

4. Quality Control

To ensure the reliability of the test results, quality control is required after each calibration, and the quality control value should be within the specified range. If it is out of the specified range, the user should check the measurement system and recalibrate it if necessary. It is recommended that users conduct regular quality

control according to their own needs to ensure the accuracy of the measurement system and the reliability of the test results.

5. Calculation

Linear calibration was used to draw the working curve. The concentration of albumin (ALB) in the sample can be calculated on the working curve based on its absorbance change value.

Reference Intervals

35.0~55.0 g/L

This reference interval is determined based on 95% distribution interval obtained from 200 healthy human specimens without related diseases, and is for reference only. It is recommended that each laboratory establish its own reference range.

Explanation of Results

If the concentration of ALB in the sample exceeds 60.00 g/L, the reduction mode of the chemistry analyzer will be used, or the high concentration sample shall be diluted with normal saline, with its report result multiplied by the dilution factor.

The professional is responsible for the review of the test result, which may be affected by the subject's age, sex, or weight. The test values within the critical range should be confirmed by remeasuring, if it is obviously beyond the reference value range or if it is still beyond the reference range after confirmation detection, the target concentration in the sample is considered to be abnormal. If the test result is inconsistent with or even contrary to the clinical situation, cause analysis should be performed correspondingly.

Limitations

1. The deviation of test results caused by interferents is less than 10% if the concentrations of the following interferents are at or below the given values:

Substances	Concentrations
Vc	0.5 g/L
Chyle	0.30%
Bilirubin	342 μ mol/L

2. The test results of this product are only for clinical reference and cannot be used as the basis for diagnosis or exclusion of cases. The clinical diagnosis and treatment of patients should be combined with their symptoms/signs, medical history, other laboratory tests and treatment response. For achieve diagnostic purposes, the test results should be combined with clinical tests, medical history and other test results.

Performance Characteristics

- The reagent blank absorbance ≤ 0.500 .
- Analytical sensitivity: at the test concentration of 40.0 g/L, the reagent absorbance change (ΔA) ≥ 0.50 .
- Accuracy: relative deviation $\leq 6.0\%$.
- Precision: within-run CV $\leq 2.0\%$, between-run relative range $\leq 5.0\%$.
- Linear Range:
[10.0, 60.0] g/L, the correlation coefficient (r) ≥ 0.990 .
[10.0, 20.0] g/L, the absolute deviation ≤ 4.0 g/L;
(20.0, 60.0] g/L, the relative deviation $\leq 10\%$.











Materials Required (but not provided)

Chemistry analyzer, Zybio Clinical Chemistry Multi-analyte Calibrator or Randox multi-analyte calibrator, Control, General lab equipment and consumable.

References

[1] Guo J, Xie J, Zhao H. Design of method comparison study and bias estimation for albumin assays[J]. Chin J Lab Med, 2000, 23:343-345.

Symbol Interpretation

	In Vitro Diagnostic Medical Device		Batch Code
	Consult Instructions for Use		Use-By Date
	Catalogue Number		Manufacturer
	Temperature Limit		Date of Manufacture
	CE marking of conformity		Authorized Representative in the European Community



Zybio Inc.

Floor 1 to Floor 5, Building 30, No. 6 of Taikang Road, Block C of Jianqiao Industrial Park, Dadukou District, 400082 Chongqing, PEOPLE'S REPUBLIC OF CHINA.

Web: www.zybio.com

E-mail: info@zybio.com

Tel: +86 (0)23 6865 5509

Fax: +86 (0)23 6869 9779



Lotus NL B.V.

Koningin Julianaplein 10, 1e Verd, 2595AA, The Hague, Netherlands.

Current Version: 02 Date of Issue: May, 2022

Instructions for Use of Alanine Aminotransferase (ALT) Kit (Enzymatic Method)

Package Specification

REF	Reagent	Systems
01.09.00.05.EC.01	R1 30 mL × 3	Zybio EXC200/220
	R2 7.5 mL × 3	
01.09.00.05.EC.03	R1 48 mL × 2	Hitachi 7180
	R2 12 mL × 2	Zybio EXC400/420

Intended Use

In vitro test for the quantitative determination of alanine aminotransferase activity in human samples (serum or plasma). Clinically, it is mainly used as an aid to diagnosis of hepatobiliary diseases.

Summary

The enzyme alanine aminotransferase (ALT) has been widely reported as present in a variety of tissues. The major source of ALT is the liver, which has led to the measurement of ALT activity for the diagnosis of hepatic diseases. Elevated serum ALT is found in hepatitis, cirrhosis, obstructive jaundice, carcinoma of the liver, and chronic alcohol abuse. ALT is only slightly elevated in patients who have an uncomplicated myocardial infarction. Although both serum aspartate aminotransferase (AST) and ALT become elevated whenever disease processes affect liver cell integrity, ALT is the more liver-specific enzyme. Moreover, elevations of ALT activity persist longer than elevations of AST activity. In patients with vitamin B6 deficiency, serum aminotransferase activity maybe decreased. The apparent reduction in aminotransferase activity may be related to decreased pyridoxal phosphate, the prosthetic group for aminotransferases, resulting in an increase in the ratio of apoenzyme to holoenzyme.

Principle

This kit uses the method recommended by the International Federation of Clinical Chemistry (IFCC):

- Alanine + α -Ketoglutaric Acid $\xrightarrow{\text{ALT}}$ Pyruvic Acid + L-Glutamic Acid
- Pyruvic Acid + NADH + H⁺ $\xrightarrow{\text{LDH}}$ L-Lactic Acid + NAD⁺ + H₂O

Oxidation of NADH to NAD⁺ causes a decrease in absorbance at 340 nm, which is directly proportional to the ALT activity in the sample.

Reagents Components and Concentration

Components	Main Constituents	Concentration
R1	Trometamol (Tris) buffer	62 mmol/L
	Nicotinamide adenine dinucleotide (NADH)	0.4 mmol/L
R2	Trometamol (Tris) buffer	512 mmol/L
	α -Ketoglutaric Acid	79.6 mmol/L
	L-Alanine	898 mmol/L
	Lactate Dehydrogenase (LDH)	≥8.5 kU/L

The components in different batches are non-interchangeable.

Storage and Validity

- The reagents should be stored at 2 - 8 °C and kept away from direct light and freezing. The unopened reagents are valid for 12 months.
- Once opened, the reagents are stable for 4 weeks at 2 - 8 °C. For reagents not in use, the cap should be tightened to avoid contamination.

Floor 1 to Floor 5, Building 30, No. 6 of Taikang Road, Block C of Jianqiao Industrial Park, Dadukou District, 400082 Chongqing, PEOPLE'S REPUBLIC OF CHINA
Tel: +86 (0)23 6865 5509 Fax: +86 (0)23 6869 9779 E-mail: info@zybio.com Web: www.zybio.com

- The production date and expiration date are available on package insert.

System Information

Hitachi 7180, Zybio EXC400/420, Zybio EXC200/220 Chemistry Analyzer. Other models can be used after verification.

Specimen Information

Non-hemolytic serum or plasma is suitable for samples, which are stable for 3 days at 2 - 8 °C. Avoid repeated freezing and thawing.

Warnings and Precautions

- For in vitro diagnostic use only. In case of contact, immediately flush the affected area with plenty of water.
- The dosage of reagents and samples can be appropriately increased or decreased according to various instruments, under the condition that the volume proportion of reagents and samples is invariable.
- The necessary precautions should be taken to use the reagents. The bottle cap should be tightened immediately after use to avoid contamination.
- When reagent becomes turbid or the blank absorbance < 1.000, the reagent is failed and should be discarded.
- All samples and reaction wastes should be treated as sources of infection. The waste liquid generated during the experiment and the used packaging materials shall be collected and disposed in accordance with relevant regulations.
- The same sample tested with reagents from different manufacturers may lead to different measured values.
- Warning: This product contains ingredients of human or animal origin. At present, there is no way to completely guarantee that it is free from infectious substances and may be contaminated during use. This product and samples shall be considered as potential sources of infection, and operators should take protective measures and comply with laboratory safety regulations. All waste shall be disposed of in accordance with local regulations.

Test Process

1. Parameters

Method	Rate Method	Sample/Reagent	6/125
Main Wavelength	340 nm	Reaction Temperature	37 °C
Sub Wavelength	405 nm	Reaction Time	10 min
Reaction Direction	-		

2. Operation

Addition	Blank	Calibration	Detection
Sample (μL)	/	/	12
Calibrator (μL)	/	12	/
Purified Water (μL)	12	/	/
Reagent 1 (μL)	200	200	200
Mix well, incubate at 37 °C for 5 min			
Reagent 2 (μL)	60	60	60
Mix well, after 2 min, accurately measure the absorbance change rate ΔA/min within 3 min.			

3. Calibration

Use Zybio Clinical Chemistry Multi-analyte Calibrator or Randox multi-analyte calibrator. Calibration cycle: Re-calibration is needed in case of replacement of the reagent batch, deviation of quality control, or the major repair and maintenance of equipment.

4. Quality Control

To ensure the reliability of the test results, quality control is required after each calibration, and the quality control value should be within the specified range. If it is out of the specified range, the user should check the measurement system and recalibrate it if necessary. It is recommended that users conduct regular quality control according to their own needs to ensure the accuracy of the measurement system and the reliability of the test results.

5. Calculation

Linear calibration was used to draw the working curve. The concentration of alanine aminotransferase (ALT) in the sample can be calculated on the working curve based on its absorbance change rate.

Reference Intervals

Male: 9~50 U/L

Female: 7~40 U/L

This reference interval is determined based on 95% distribution interval obtained from 200 healthy males and 200 healthy females specimens without related diseases, and is for reference only. It is recommended that each laboratory establish its own reference range.

Explanation of Results

If the concentration of ALT in the sample exceeds 1000 U/L, the reduction mode of the chemistry analyzer will be used, or the high concentration sample shall be diluted with normal saline, with its report result multiplied by the dilution factor.

The professional is responsible for the review of the test result, which may be affected by the subject's age, sex, or weight. The test values within the critical range should be confirmed by remeasuring, if it is obviously beyond the reference value range or if it is still beyond the reference range after confirmation detection, the target concentration in the sample is considered to be abnormal. If the test result is inconsistent with or even contrary to the clinical situation, cause analysis should be performed correspondingly.

Limitations

1. The deviation of test results caused by interferents is less than 10% if the concentrations of the following interferents are at or below the given values:

Substances	Concentrations
Hemoglobin	5 g/L
Chyle	0.30%
Bilirubin	300 μ mol/L
Triglyceride	11.3 mmol/L

2. The test results of this product are only for clinical reference and cannot be used as the basis for diagnosis or exclusion of cases. The clinical diagnosis and treatment of patients should be combined with their symptoms/signs, medical history, other laboratory tests and treatment response. For achieve diagnostic purposes, the test results should be combined with clinical tests, medical history and other test results.

Performance Characteristics

- The reagent blank absorbance ≥ 1.000 ; the reagent blank absorbance change rate ($\Delta A/\text{min}$) ≤ 0.004 .
- Analytical sensitivity: at the test concentration of 130 U/L, the reagent absorbance change rate ($\Delta A/\text{min}$) ≥ 0.01 .
- Accuracy: relative deviation $\leq 10\%$.

4. Precision: within-run CV $\leq 5\%$, between-run relative range $\leq 10\%$.

5. Linear Range:

[5, 1000] U/L, the correlation coefficient (r) ≥ 0.990 .

[5, 40] U/L, the absolute deviation ≤ 4 U/L;

(40, 1000] U/L, the relative deviation $\leq 10\%$.

Materials Required (but not provided)

Chemistry analyzer, Zybion Clinical Chemistry Multi-analyte Calibrator or Randox multi-analyte calibrator, Control, General lab equipment and consumable.

References

[1] Prati D, Taioli E, Zanella A, et al. Updated definitions of healthy ranges for serum alanine aminotransferase levels[J]. Ann Intern Med, 2002, 137:1-10.

Symbol Interpretation

IVD	In Vitro Diagnostic Medical Device	LOT	Batch Code
	Consult Instructions for Use		Use-By Date
REF	Catalogue Number		Manufacturer
	Temperature Limit		Date of Manufacture
CE	CE marking of conformity	EC REP	Authorized Representative in the European Community



Zybio Inc.

Floor 1 to Floor 5, Building 30, No. 6 of Taikang Road, Block C of Jianqiao Industrial Park, Dadukou District, 400082 Chongqing, PEOPLE'S REPUBLIC OF CHINA.

Web: www.zybio.com

E-mail: info@zybio.com

Tel: +86 (0)23 6865 5509

Fax: +86 (0)23 6869 9779



Lotus NL B.V.

Koningin Julianaplein 10, 1e Verd, 2595AA, The Hague, Netherlands.

Current Version: 03

Date of Issue: April, 2023

Instructions for Use of Aspartate Aminotransferase (AST) Kit (Enzymatic Method)

Package Specification

REF	Reagent	Systems
01.09.00.16.EC.01	R1 30 mL × 3	Zybio EXC200/220
	R2 7.5 mL × 3	
01.09.00.16.EC.02	R1 48 mL × 2	Hitachi 7180
	R2 12 mL × 2	Zybio EXC400/420

Intended Use

In vitro test for the quantitative determination of aspartate aminotransferase activity in human samples (serum or plasma). Clinically, it is mainly used as an aid to diagnosis of viral hepatitis, obstructive jaundice, and myocardial infarction.

Summary

The enzyme aspartate aminotransferase (AST) is widely distributed in tissue, principally hepatic, cardiac, muscle, and kidney. Elevated serum levels are found in diseases involving these tissues. Hepatobiliary diseases, such as cirrhosis, metastatic carcinoma, and viral hepatitis also increase serum AST levels. Following myocardial infarction, serum AST is elevated and reaches a peak two days after onset. In patients undergoing renal dialysis or those with vitamin B6 deficiency, serum AST may be decreased. The apparent reduction in AST may be related to decreased pyridoxal phosphate, the prosthetic group for AST, resulting in an increase in the ratio of apoenzyme to holoenzyme. Two isoenzymes of AST have been detected, cytoplasmic and mitochondrial. Only the cytoplasmic isoenzyme occurs in normal serum, while the mitochondrial, together with the cytoplasmic isoenzyme, has been detected in the serum of patients with coronary and hepatobiliary disease.

Principle

This kit uses the method recommended by the International Federation of Clinical Chemistry (IFCC):

- Aspartic Acid + α -Ketoglutaric Acid $\xrightarrow{\text{AST}}$ Oxaloacetic Acid + L-Glutamic Acid
- Oxaloacetic Acid + NADH + H⁺ $\xrightarrow{\text{MDH}}$ L-Lactic Acid + NAD⁺ + H₂O

Oxidation of NADH to NAD⁺ causes a decrease in absorbance at 340 nm, which is directly proportional to the AST activity in the sample.

Reagents Components and Concentration

Components	Main Constituents	Concentration
R1	Trometamol (Tris) buffer	62 mmol/L
	Nicotinamide adenine dinucleotide (NADH)	0.4 mmol/L
R2	Trometamol (Tris) buffer	439 mmol/L
	α -Ketoglutaric Acid	37.1 mmol/L
	L-Aspartic Acid	>800 mmol/L
	Malate Dehydrogenase (MDH)	>2.5 kU/L

The components in different batches are non-interchangeable.

Storage and Validity

- The reagents should be stored at 2 - 8 °C and kept away from direct light and freezing. The unopened reagents are valid for 12 months.
- Once opened, the reagents are stable for 4 weeks at 2 - 8 °C. For reagents not in use, the cap should be tightened to avoid contamination.

Floor 1 to Floor 5, Building 30, No. 6 of Taikang Road, Block C of Jianqiao Industrial Park, Dadukou District, 400082 Chongqing, PEOPLE'S REPUBLIC OF CHINA
Tel: +86 (0)23 6865 5509 Fax: +86 (0)23 6869 9779 E-mail: info@zybio.com Web: www.zybio.com

- The production date and expiration date are available on package insert.

System Information

Hitachi 7180, Zybio EXC400/420, Zybio EXC200/220 Chemistry Analyzer. Other models can be used after verification.

Specimen Information

Non-hemolytic serum or plasma is suitable for samples, which are stable for 3 days at 2 - 8 °C. Avoid repeated freezing and thawing.

Warnings and Precautions

- For in vitro diagnostic use only. In case of contact, immediately flush the affected area with plenty of water.
- The dosage of reagents and samples can be appropriately increased or decreased according to various instruments, under the condition that the volume proportion of reagents and samples is invariable.
- The necessary precautions should be taken to use the reagents. The bottle cap should be tightened immediately after use to avoid contamination.
- When reagent becomes turbid or the blank absorbance < 1.000, the reagent is failed and should be discarded.
- All samples and reaction wastes should be treated as sources of infection. The waste liquid generated during the experiment and the used packaging materials shall be collected and disposed in accordance with relevant regulations.
- The same sample tested with reagents from different manufacturers may lead to different measured values.
- Warning: This product contains ingredients of human or animal origin. At present, there is no way to completely guarantee that it is free from infectious substances and may be contaminated during use. This product and samples shall be considered as potential sources of infection, and operators should take protective measures and comply with laboratory safety regulations. All waste shall be disposed of in accordance with local regulations.

Test Process

1. Parameters

Method	Rate Method	Sample/Reagent	6/125
Main Wavelength	340 nm	Reaction Temperature	37 °C
Sub Wavelength	405 nm	Reaction Time	10 min
Reaction Direction	-		

2. Operation

Addition	Blank	Calibration	Detection
Sample (μL)	/	/	12
Calibrator (μL)	/	12	/
Purified Water (μL)	12	/	/
Reagent 1 (μL)	200	200	200
Mix well, incubate at 37 °C for 5 min			
Reagent 2 (μL)	50	50	50
Mix well, after 2 min, measure the average absorbance change rate $\Delta A/\text{min}$ within 3 min.			

3. Calibration

Use Zybio Clinical Chemistry Multi-analyte Calibrator or Randox multi-analyte calibrator. Calibration cycle: Re-calibration is needed in case of replacement of the reagent batch, deviation of quality control, or the major repair and maintenance of equipment.

4. Quality Control

To ensure the reliability of the test results, quality control is required after each calibration, and the quality control value should be within the specified range. If it is out of the specified range, the user should check the measurement system and recalibrate it if necessary. It is recommended that users conduct regular quality control according to their own needs to ensure the accuracy of the measurement system and the reliability of the test results.

5. Calculation

Linear calibration was used to draw the working curve. The concentration of aspartate aminotransferase (AST) in the sample can be calculated on the working curve based on its absorbance change rate.

Reference Intervals

≤ 40 U/L

This reference interval is determined based on 95% distribution interval obtained from 200 healthy human specimens without related diseases, and is for reference only. It is recommended that each laboratory establish its own reference range.

Explanation of Results

If the concentration of AST in the sample exceeds 1000 U/L, the reduction mode of the chemistry analyzer will be used, or the high concentration sample shall be diluted with normal saline, with its report result multiplied by the dilution factor.

The professional is responsible for the review of the test result, which may be affected by the subject's age, sex, or weight. The test values within the critical range should be confirmed by remeasuring, if it is obviously beyond the reference value range or if it is still beyond the reference range after confirmation detection, the target concentration in the sample is considered to be abnormal. If the test result is inconsistent with or even contrary to the clinical situation, cause analysis should be performed correspondingly.

Limitations

1. The deviation of test results caused by interferents is less than 10% if the concentrations of the following interferents are at or below the given values:

Substances	Concentrations
Chyle	0.30%
Bilirubin	300 μmol/L

2. The test results of this product are only for clinical reference and cannot be used as the basis for diagnosis or exclusion of cases. The clinical diagnosis and treatment of patients should be combined with their symptoms/signs, medical history, other laboratory tests and treatment response. For achieve diagnostic purposes, the test results should be combined with clinical tests, medical history and other test results.

Performance Characteristics

- The reagent blank absorbance ≥ 1.000; the reagent blank absorbance change rate ($\Delta A/\text{min}$) ≤ 0.004.
- Analytical sensitivity: at the test concentration of 130.0 U/L, the reagent absorbance change rate ($\Delta A/\text{min}$) ≥ 0.01.
- Accuracy: relative deviation ≤ 10%.
- Precision: within-run CV ≤ 5%, between-run relative range ≤ 10%.
- Linear Range:
 - [10, 1000] U/L, the correlation coefficient (r) ≥ 0.990.
 - [10, 100] U/L, the absolute deviation ≤ 10 U/L;
 - (100, 1000] U/L, the relative deviation ≤ 10%.

Materials Required (but not provided)

Chemistry analyzer, Zybio Clinical Chemistry Multi-analyte Calibrator or Randox multi-analyte calibrator, Control, General lab equipment and consumable.

References

- [1] Abdalla D. Clinical chemistry: theory, analysis, correlations[J]. Revista Brasileira de Ciências Farmacêuticas, 2003, 39:348-349.
- [2] Tietz N. Fundamentals of clinical chemistry[M]. Saunders, 1987.

Symbol Interpretation

IVD	In Vitro Diagnostic Medical Device	LOT	Batch Code
	Consult Instructions for Use		Use-By Date
REF	Catalogue Number		Manufacturer
	Temperature Limit		Date of Manufacture
CE	CE marking of conformity	EC REP	Authorized Representative in the European Community



Zybio Inc.
 Floor 1 to Floor 5, Building 30, No. 6 of Taikang Road, Block C of Jianqiao Industrial Park, Dadukou District, 400082 Chongqing, PEOPLE'S REPUBLIC OF CHINA.
 Web: www.zybio.com
 E-mail: info@zybio.com
 Tel: +86 (0)23 6865 5509
 Fax: +86 (0)23 6869 9779



Lotus NL B.V.
 Koningin Julianaplein 10, 1e Verd, 2595AA, The Hague, Netherlands.

Current Version: 02 Date of Issue: May, 2022

Instructions for Use of Calcium (Ca) Kit (Arsenazo III Method)

Package Specification

REF	Reagent	Systems
01.09.0C.01.EC.01	R 30 mL x 6	Zybio EXC200/220
01.09.0C.01.EC.02	R 60 mL x 2	Hitachi 7180 Zybio EXC400/420

Intended Use

In vitro test for the quantitative determination of calcium (Ca) concentration in human samples (serum or plasma). Clinically, it is mainly used as an aid to diagnosis of calcium metabolism disorders.

Summary

Calcium is the most abundant mineral element in the body with about 99% in the bones primarily as hydroxyapatite. The remaining calcium is distributed between the various tissues and the extracellular fluids where it performs a vital role for many life sustaining processes. Among the extra skeletal functions of calcium are involvement in blood coagulation, neuromuscular conduction, excitability of skeletal and cardiac muscle, enzyme activation, and the preservation of cell membrane integrity and permeability. Serum calcium levels and hence the body content are controlled by parathyroid hormone (PTH), calcitonin, and vitamin D. An imbalance in any of these modulators leads to alterations of the body and serum calcium levels. Increases in serum PTH or vitamin D are usually associated with hypercalcemia. Increased serum calcium levels may also be observed in multiple myeloma and other neoplastic diseases. Hypocalcemia may be observed e.g. in hypoparathyroidism, nephrosis, and pancreatitis.

Principle

The Arsenazo III is combined with calcium ions, forming a purple-colored complex. The color of the complex is proportional to the concentration of calcium ion in the sample, which can be calculated by measuring the absorbance change at 660 nm.

Reagents Components and Concentration

Components	Main Constituents	Concentration
R	Arsenazo III	129 µmol/L
	MES Buffer	4.25 g/L
	Surfactant	0.2% (v/v)

The components in different batches are non-interchangeable.

Storage and Validity

- The reagents should be stored at 2 - 8 °C and kept away from direct light and freezing. The unopened reagents are valid for 12 months.
- Once opened, the reagents are stable for 4 weeks at 2 - 8 °C. For reagents not in use, the cap should be tightened to avoid contamination.
- The production date and expiration date are available on package insert.

System Information

Hitachi 7180, Zybio EXC400/420, Zybio EXC200/220 Chemistry Analyzer. Other models can be used after verification.

Specimen Information

- Fresh and nonhemolytic serum or plasma (heparin) is suitable for samples.
- Samples should be analyzed as soon as possible after collection, which can be

Floor 1 to Floor 5, Building 30, No. 6 of Taikang Road, Block C of Jianqiao Industrial Park, Dadukou District, 400082 Chongqing, PEOPLE'S REPUBLIC OF CHINA
Tel: +86 (0)23 6865 5509 Fax: +86 (0)23 6869 9779 E-mail: info@zybio.com

stable for 2 days at 20 - 25 °C, for 14 days at 2 - 8 °C, and for 3 months at - 20 °C. Repeated freezing and thawing should be avoided.

Warnings and Precautions

- For in vitro diagnostic use only. In case of contact, immediately flush the affected area with plenty of water.
- The dosage of reagents and samples can be appropriately increased or decreased according to various instruments, under the condition that the volume proportion of reagents and samples is invariable.
- Strict measures shall be taken to avoid contamination since calcium ion is almost omnipresent.
- When reagent becomes turbid or the blank absorbance > 1.500, the reagent is failed and should be discarded.
- All samples and reaction wastes should be treated as sources of infection. The waste liquid generated during the experiment and the used packaging materials shall be collected and disposed in accordance with relevant regulations.
- Trace chelating agents (such as EDTA) present in the detergent can hinder the generation of chromogens. It is recommended to use disposable tubes and pipettes, etc.
- The same sample tested with reagents from different manufacturers may lead to different measured values.
- Warning: This product contains ingredients of human or animal origin. At present, there is no way to completely guarantee that it is free from infectious substances and may be contaminated during use. This product and samples shall be considered as potential sources of infection, and operators should take protective measures and comply with laboratory safety regulations. All waste shall be disposed of in accordance with local regulations.

Test Process

1. Parameters

Method	End-Point Method	Sample/Reagent	1/100
Main Wavelength	660 nm	Reaction Temperature	37 °C
Sub Wavelength	700 nm	Reaction Direction	+

2. Operation

Addition	Blank	Calibration	Detection
Sample (µL)	/	/	3
Calibrator (µL)	/	3	/
Purified Water (µL)	3	/	/
Reagent (µL)	300	300	300
Mix well, incubate at 37 °C for 2 min, then zero the system at 660 nm as blank and measure absorbance A.			

3. Calibration

Use Zybio Clinical Chemistry Multi-analyte Calibrator or Randox multi-analyte calibrator. Calibration cycle: Re-calibration is needed in case of replacement of the reagent batch, deviation of quality control, or the major repair and maintenance of equipment.

4. Quality Control

To ensure the reliability of the test results, quality control is required after each calibration, and the quality control value should be within the specified range. If it is out of the specified range, the user should check the measurement system and recalibrate it if necessary. It is recommended that users conduct regular quality control according to their own needs to ensure the accuracy of the measurement system and the reliability of the test results.

Dadukou District, 400082 Chongqing, PEOPLE'S REPUBLIC OF CHINA
Web: www.zybio.com

5. Calculation

Linear calibration was used to draw the working curve. The concentration of calcium ion (Ca) in the sample can be calculated on the working curve based on its absorbance change value.

Reference Intervals

Adults Serum: 2.10 - 2.60 mmol/L

Children Serum: 2.50 - 3.00 mmol/L

This reference interval is determined based on 95% distribution interval obtained from 210 healthy human specimens without related diseases per group, and is for reference only. It is recommended that each laboratory establish its own reference range.

Explanation of Results

If the concentration of Ca in the sample exceeds 4.00 mmol/L, the reduction mode of the chemistry analyzer will be used, or the high concentration sample shall be diluted with normal saline, with its report result multiplied by the dilution factor.

The professional is responsible for the review of the test result, which may be affected by the subject's age, sex, or weight. The test values within the critical range should be confirmed by remeasuring, if it is obviously beyond the reference value range or if it is still beyond the reference range after confirmation detection, the target concentration in the sample is considered to be abnormal. If the test result is inconsistent with or even contrary to the clinical situation, cause analysis should be performed correspondingly.

Limitations

1. The deviation of test results caused by interferences is within $\pm 10\%$ if the concentrations of the following interferences are at or below the given values:

Substances	Concentrations
Bilirubin	280 $\mu\text{mol/L}$
Mg ²⁺	3 mmol/L
K ⁺	8 mmol/L
Na ⁺	180 mmol/L

2. The test results of this product are only for clinical reference and cannot be used as the basis for diagnosis or exclusion of cases. The clinical diagnosis and treatment of patients should be combined with their symptoms/signs, medical history, other laboratory tests and treatment response. For achieve diagnostic purposes, the test results should be combined with clinical tests, medical history and other test results.

Performance Characteristics

- The reagent blank absorbance ≤ 1.500 .
- Analytical sensitivity: at the test concentration of 2.50 mmol/L, the absorbance change (ΔA) ≥ 0.20 .
- Accuracy: relative deviation $\leq 5\%$.
- Precision: within-run CV $\leq 3\%$, between-run relative range $\leq 5\%$.
- Linear range:
[1.00, 4.00] mmol/L, the correlation coefficient (r) ≥ 0.990 .
Within the specified test range, the relative deviation $\leq 10\%$.

Materials Required (but not provided)

Chemistry analyzer, Zybio Clinical Chemistry Multi-analyte Calibrator or Randox multi-analyte calibrator, Control, General lab equipment and consumable.

References

[1] Massry S, Coburn J, Chapman L, et al. Role of serum Ca, parathyroid hormone, and NaCl infusion on renal Ca and Na clearances[J]. Am J Physiol, 1968, 214:1403-1409.

Symbol Interpretation

IVD	In Vitro Diagnostic Medical Device	LOT	Batch Code
	Consult Instructions for Use		Use-By Date
REF	Catalogue Number		Manufacturer
	Temperature Limit		Date of Manufacture
CE	CE marking of conformity	EC REP	Authorized Representative in the European Community



Zybio Inc.

Floor 1 to Floor 5, Building 30, No. 6 of Taikang Road, Block C of Jianqiao Industrial Park, Dadukou District, 400082 Chongqing, PEOPLE'S REPUBLIC OF CHINA.

Web: www.zybio.com

E-mail: info@zybio.com

Tel: +86 (0)23 6865 5509

Fax: +86 (0)23 6869 9779



Lotus NL B.V.

Koningin Julianaplein 10, 1e Verd, 2595AA, The Hague, Netherlands.

Current Version: 02 Date of Issue: May, 2022

Instructions for Use of Creatine Kinase (CK) Kit (Rate Method)

Package Specification

REF	Reagent	Systems
01.09.04.03.EC.01	R1 30 mL x 3 R2 7.5 mL x 3	Zybio EXC200/220
01.09.04.03.EC.02	R1 48 mL x 2 R2 12 mL x 2	Hitachi 7180 Zybio EXC400/420

Intended Use

In vitro test for the quantitative determination of the catalytic activity concentration of creatine kinase (CK) in human samples (serum or plasma). Clinically, serum creatine kinase levels increased in certain tissue damage or diseases, such as myocardial infarction, muscular dystrophy, acute cerebrovascular accident, etc.

Summary

Creatine kinase (CK), also known as creatine phosphokinase. The contents of creatine kinase in skeletal muscle, myocardium, and smooth muscle were more, followed by brain tissue, and the contents in gastrointestinal tract, lung and kidney were less. CK is an important kinase that is directly related to intracellular energy transport, muscle contraction, and ATP regeneration. When the striated muscle of the human body is damaged and necrotic, creatine kinase is released into the blood and abnormally elevated during detection. However, the increase of creatine kinase lacks specificity.

If it is combined with the increase of creatine kinase isoenzyme and troponin, it shall be judged that the myocardium is damaged. If abnormal elevation of myoglobin is also combined, it shall be judged that skeletal muscle injury occurs.

Principle

This kit uses rate method (improved based on the method recommended by IFCC) to determine the catalytic activity concentration of CK in samples. Creatine kinase (CK) catalyzes the conversion of phosphocreatine to creatine, while ADP is phosphorylated to ATP. Hexokinase (HK) catalyzes the reaction between ATP and glucose to generate glucose-6-phosphate. Glucose-6-phosphate dehydrogenase (G6PDH) catalyzes glucose-6-phosphate to generate 6-phosphogluconic acid, and simultaneously converts oxidative nicotinamide adenine dinucleotide phosphate (NADP⁺) into reduced nicotinamide adenine dinucleotide phosphate (NADPH), causing the increase of the absorbance. The rate of increase is directly proportional to the catalytic activity concentration of CK in the sample. By continuously monitor the absorbance change rate, the catalytic activity concentration of CK in the sample can be calculated from the calibration curve generated by the calibrator treated in the same manner.

1. Phosphocreatine + ADP $\xrightarrow{\text{CK}}$ Creatine + ATP
2. ATP + Glucose $\xrightarrow{\text{HK}}$ Glucose-6-phosphate + ADP
3. Glucose-6-phosphate + NADP⁺ $\xrightarrow{\text{G6PDH}}$ 6-phosphogluconic acid + NADPH + H⁺

Reagents Components and Concentration

Components	Main Constituents	Concentration
R1	D-Glucose	20-30 mmol/L
	Nicotinamide adenine dinucleotide phosphate oxidized form (NADP ⁺)	1.5-2.0 g/L
	Hexokinase	4-8 kU/L
	Imidazole buffer	100 mmol/L

R2	Glucose-6-phosphate dehydrogenase (G6PDH)	12-16 kU/L
	Phosphocreatine	80-120 mmol/L
	Adenosine-5'-diphosphate (ADP) potassium salt	7-9 mmol/L

The components in different batches are non-interchangeable.

Storage and Validity

1. The reagents should be stored at 2 - 8 °C and kept away from freezing. The unopened reagents are valid for 18 months.
2. Once opened, the reagents are stable for 30 days at 2 - 8 °C. For reagents not in use, the cap should be tightened to avoid contamination.
3. The production date and expiration date are available on package insert.

System Information

Hitachi 7180, Zybio EXC400/420, Zybio EXC200/220 Chemistry Analyzer. Other models can be used after verification.

Specimen Information

Serum or plasma (heparin for anticoagulation) is suitable for samples, which shall be separated in time after collection to avoid hemolysis. Samples containing EDTA, citrate, and chloride should not be used. Samples are stable for 1 day at 2 - 8 °C and 30 days at - 20 °C. Avoid repeated freezing and thawing.

Warnings and Precautions

1. For in vitro diagnostic use only. In case of contact, immediately flush the affected area with plenty of water.
2. The dosage of reagents and samples can be appropriately increased or decreased according to various instruments, under the condition that the volume proportion of reagents and samples is invariable.
3. The necessary precautions should be taken to use the reagents. The bottle cap should be tightened immediately after use to avoid contamination.
4. When the blank absorbance > 0.400, the reagent is failed and should be discarded.
5. All samples and reaction wastes should be treated as sources of infection. The waste liquid generated during the experiment and the used packaging materials shall be collected and disposed in accordance with relevant regulations.
6. The same sample tested with reagents from different manufacturers may lead to different measured values.
7. Warning: This product contains ingredients of human or animal origin. At present, there is no way to completely guarantee that it is free from infectious substances and may be contaminated during use. This product and samples shall be considered as potential sources of infection, and operators should take protective measures and comply with laboratory safety regulations. All waste shall be disposed of in accordance with local regulations.

Test Process

1. Parameters

Method	Rate Method	Sample/Reagent	1/25
Main Wavelength	340 nm	Reaction Temperature	37 °C
Sub Wavelength	546 nm	Reaction Time	10 min
Reaction Direction	+		

2. Operation

Addition	Blank	Calibration	Detection
Sample (μL)	/	/	10
Calibrator (μL)	/	10	/
Purified Water (μL)	10	/	/
Reagent 1 (μL)	200	200	200
Mix well, incubate at 37 °C for 5 min			
Reagent 2 (μL)	50	50	50
Mix well, after adding R2 for 45 s, continuously monitor the absorbance change within 4 min and 15 s, and calculate the absorbance change rate ΔA/min.			

3. Calibration

Use Zybio Clinical Chemistry Multi-analyte Calibrator or Randox multi-analyte calibrator. Calibration cycle: Re-calibration is needed in case of replacement of the reagent batch, deviation of quality control, or the major repair and maintenance of equipment.

4. Quality Control

To ensure the reliability of the test results, quality control is required after each calibration, and the quality control value should be within the specified range. If it is out of the specified range, the user should check the measurement system and recalibrate it if necessary. It is recommended that users conduct regular quality control according to their own needs to ensure the accuracy of the measurement system and the reliability of the test results.

5. Calculation

Linear calibration was used to draw the working curve. The catalytic activity concentration of creatine kinase (CK) in the sample can be calculated on the working curve based on its absorbance change rate.

Reference Intervals

Male: 38~174 U/L Female: 26~140 U/L

This reference interval is determined based on 95% distribution interval obtained from 200 healthy human specimens without related diseases per group, and is for reference only. It is recommended that each laboratory establish its own reference range.

Explanation of Results

If the catalytic activity concentration of CK in the sample exceeds 1000 U/L, the reduction mode of the chemistry analyzer will be used, or the high concentration sample shall be diluted with normal saline, with its report result multiplied by the dilution factor.

The professional is responsible for the review of the test result, which may be affected by the subject's age, sex, or weight. The test values within the critical range should be confirmed by remeasuring, if it is obviously beyond the reference value range or if it is still beyond the reference range after confirmation detection, the target concentration in the sample is considered to be abnormal. If the test result is inconsistent with or even contrary to the clinical situation, cause analysis should be performed correspondingly.

Limitations

1. The deviation of test results caused by interferents is ≤ 10% if the concentrations of the following interferents are at or below the given values:

Substances	Concentrations
Vc	1.0 g/L
Hemoglobin	5 g/L

Bilirubin	342 μmol/L
Triglyceride	10 mmol/L

2. The test results of this product are only for clinical reference and cannot be used as the basis for diagnosis or exclusion of cases. The clinical diagnosis and treatment of patients should be combined with their symptoms/signs, medical history, other laboratory tests and treatment response. For achieve diagnostic purposes, the test results should be combined with clinical tests, medical history and other test results.

Performance Characteristics

- The reagent blank absorbance ≤ 0.400; the reagent blank absorbance change rate (ΔA/min) ≤ 0.002.
- Analytical sensitivity: at the test catalytic activity concentration of 100 U/L, the reagent absorbance change rate (ΔA/min) > 0.002.
- Accuracy: relative deviation ≤ 10%.
- Precision: within-run CV ≤ 5%, between-run relative range ≤ 10%.
- Linear Range:
[25, 1000] U/L, the correlation coefficient (r) ≥ 0.990.
[25, 100] U/L, the absolute deviation ≤ 10 U/L;
[100, 1000] U/L, the relative deviation ≤ 10%.

Materials Required (but not provided)

Chemistry analyzer, Zybio Clinical Chemistry Multi-analyte Calibrator or Randox multi-analyte calibrator, Control, General lab equipment and consumable.

References

[1] Kitzenberg D, Colgan S, Glover L. Creatine kinase in ischemic and inflammatory disorders[J]. Clin Transl Med, 2016, 5:31.

Symbol Interpretation

IVD	In Vitro Diagnostic Medical Device	LOT	Batch Code
	Consult Instructions for Use		Use-By Date
REF	Catalogue Number		Manufacturer
	Temperature Limit		Date of Manufacture
CE	CE marking of conformity	EC REP	Authorized Representative in the European Community



Zybio Inc.
Floor 1 to Floor 5, Building 30, No. 6 of Taikang Road, Block C of Jianqiao Industrial Park, Dadukou District, 400082 Chongqing, PEOPLE'S REPUBLIC OF CHINA.
Web: www.zybio.com
E-mail: info@zybio.com
Tel: +86 (0)23 6865 5509
Fax: +86 (0)23 6869 9779



Lotus NL B.V.
Koningin Julianaplein 10, 1e Verd, 2595AA, The Hague, Netherlands.

Current Version: 02 Date of Issue: May, 2022

Instructions for Use of Gamma-Glutamyl Transferase (GGT) Kit (Enzymatic Method)

Package Specification

REF	Reagent	Systems
01.09.00.03.EC.01	R1 30 mL × 3	Zybio EXC200/220
	R2 7.5 mL × 3	
01.09.00.03.EC.02	R1 48 mL × 2	Hitachi 7180
	R2 12 mL × 2	Zybio EXC400/420

Intended Use


In vitro test for the quantitative determination of the catalytic activity concentration of γ-glutamyl transferase in human samples (serum or plasma). Clinically, it is mainly used as an aid to diagnosis of hepatobiliary diseases.

Summary

γ-Glutamyl transferase is used in the diagnosis and monitoring of hepatobiliary diseases. Enzymatic activity of GGT is often the only parameter with increased values when testing for such diseases, and is one of the most sensitive indicators known. γ-Glutamyl transferase is also a sensitive screening test for occult alcoholism. Elevated GGT activities are found in the serum of patients requiring long-term medication with phenobarbital and phenytoin. In 1969, Szasz published the first kinetic procedure for GGT in serum using γ-glutamyl-p-nitroanilide as substrate and glycylglycine as acceptor. In order to circumvent the poor solubility of γ-glutamyl-p-nitroanilide, Persijn and van der Slik investigated various derivatives and found the water soluble substrate L-γ-glutamyl-3-carboxy-4-nitroanilide to be superior in terms of stability and solubility. The results correlate with those derived using the original substrate. In 2002, the International Federation of Clinical Chemistry (IFCC) recommended the standardized method for determining GGT including optimization of substrate concentrations, employment of NaOH, glycylglycine buffer and sample start.

Principle

The kit uses a modified version of the method recommended by the International Federation of Clinical Chemistry (IFCC):

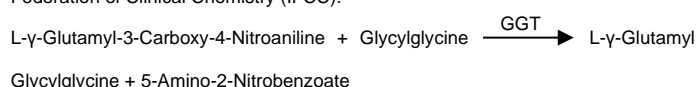

This causes an increase in absorbance at 405 nm, which is directly proportional to the catalytic activity concentration of GGT in the sample.

Reagents Components and Concentration

Components	Main Constituents	Concentration
R1	Glycylglycine	127.8 mmol/L
	Trometamol (Tris) buffer	154.6 mmol /L
R2	L-γ-Glutamyl-3-Carboxy-4-Nitroaniline	6 g/L

The components in different batches are non-interchangeable.

Storage and Validity

- The reagents should be stored at 2 - 8 °C and kept away from direct light and freezing. The unopened reagents are valid for 12 months.
- Once opened, the reagents are stable for 4 weeks at 2 - 8 °C. For reagents not in use, the cap should be tightened to avoid contamination.
- The production date and expiration date are available on package insert.

System Information

Hitachi 7180, Zybio EXC400/420, Zybio EXC200/220 Chemistry Analyzer. Other models can be used after verification.

Specimen Information

Non-hemolytic serum or plasma (EDTA for anticoagulation) is suitable for samples. The γ-glutamyl transferase in samples is stable for 7 days at 2 - 8 °C.

Warnings and Precautions

- For in vitro diagnostic use only. In case of contact, immediately flush the affected area with plenty of water.
- The dosage of reagents and samples can be appropriately increased or decreased according to various instruments, under the condition that the volume proportion of reagents and samples is invariable.
- The necessary precautions should be taken to use the reagents. The bottle cap should be tightened immediately after use to avoid contamination.
- When reagent becomes turbid or the blank absorbance > 0.800, the reagent is failed and should be discarded.
- All samples and reaction wastes should be treated as sources of infection. The waste liquid generated during the experiment and the used packaging materials shall be collected and disposed in accordance with relevant regulations.
- Considering the reaction solution turbidity caused by heparin and inhibition of GGT by citrate, oxalate, and fluoride, plasma with these substances as anticoagulant is not suitable for GGT determination.
- The same sample tested with reagents from different manufacturers may lead to different measured values.
- Warning: This product contains ingredients of human or animal origin. At present, there is no way to completely guarantee that it is free from infectious substances and may be contaminated during use. This product and samples shall be considered as potential sources of infection, and operators should take protective measures and comply with laboratory safety regulations. All waste shall be disposed of in accordance with local regulations.

Test Process

1. Parameters

Method	Rate Method	Sample/Reagent	1/10
Main Wavelength	405 nm	Reaction Temperature	37 °C
Sub Wavelength	505 nm	Reaction Time	10 min
Reaction Direction	+		

2. Operation

Addition	Blank	Calibration	Detection
Sample (μL)	/	/	25
Calibrator (μL)	/	25	/
Purified Water (μL)	25	/	/
Reagent 1 (μL)	200	200	200
Mix well, incubate at 37 °C for 3 ~ 5 min			
Reagent 2 (μL)	50	50	50
After 1 min, continuously monitor the absorbance change within 2 min, and calculate the absorbance change rate ΔA/min.			

3. Calibration

Use Zybio Clinical Chemistry Multi-analyte Calibrator or Randox multi-analyte calibrator. Calibration cycle: Re-calibration is needed in case of replacement of the reagent batch, deviation of quality control, or the major repair and maintenance of equipment.

4. Quality Control

To ensure the reliability of the test results, quality control is required after each calibration, and the quality control value should be within the specified range. If it is out of the specified range, the user should check the measurement system and recalibrate it if necessary. It is recommended that users conduct regular quality control according to their own needs to ensure the accuracy of the measurement system and the reliability of the test results.

5. Calculation

Linear calibration was used to draw the working curve. The catalytic activity concentration of γ -glutamyl transferase (GGT) in the sample can be calculated on the working curve based on its absorbance change rate.

Reference Intervals

Male: 11 - 50 U/L Female: 7 - 32 U/L

This reference interval is determined based on 95% distribution interval obtained from 200 healthy human specimens without related diseases per group, and is for reference only. It is recommended that each laboratory establish its own reference range.

Explanation of Results

If the catalytic activity concentration of GGT in the sample exceeds 600 U/L, the reduction mode of the chemistry analyzer will be used, or the high concentration sample shall be diluted with normal saline, with its report result multiplied by the dilution factor.

The professional is responsible for the review of the test result, which may be affected by the subject's age, sex, or weight. The test values within the critical range should be confirmed by remeasuring, if it is obviously beyond the reference value range or if it is still beyond the reference range after confirmation detection, the target concentration in the sample is considered to be abnormal. If the test result is inconsistent with or even contrary to the clinical situation, cause analysis should be performed correspondingly.

Limitations

1. The deviation of test results caused by interferents is less than 10% if the concentrations of the following interferents are at or below the given values:

Substances	Concentrations
Hemoglobin	5 g/L
Bilirubin	684 μ mol/L
Triglyceride	10 g/L

2. The test results of this product are only for clinical reference and cannot be used as the basis for diagnosis or exclusion of cases. The clinical diagnosis and treatment of patients should be combined with their symptoms/signs, medical history, other laboratory tests and treatment response. For achieve diagnostic purposes, the test results should be combined with clinical tests, medical history and other test results.

Performance Characteristics

- The reagent blank absorbance ≤ 0.800 , the reagent blank absorbance change rate ($\Delta A/\text{min}$) ≤ 0.005 .
- Analytical sensitivity: at the test catalytic activity concentration of 50 U/L, the absorbance change rate ($\Delta A/\text{min}$) ≥ 0.010 .
- Accuracy: relative deviation $\leq 10\%$.
- Precision: within-run CV $\leq 5\%$, between-run relative range $\leq 10\%$.

5. Linear range:

[10, 600] U/L, the correlation coefficient (r) ≥ 0.990 .

[10, 50] U/L, the absolute deviation ≤ 5 U/L;

(50, 600] U/L, the relative deviation $\leq 10\%$.

Materials Required (but not provided)

Chemistry analyzer, Zybion Clinical Chemistry Multi-analyte Calibrator or Randox multi-analyte calibrator, Control, General lab equipment and consumable.

References

- Szasz G. A kinetic photometric method for serum gamma-glutamyl transpeptidase[J]. Clin Chem, 1969, 15:124-136.
- Schumann G, Bonora R, Ceriotti F, et al. IFCC primary reference procedures for the measurement of catalytic activity concentrations of enzymes at 37 degrees C. International Federation of Clinical Chemistry and Laboratory Medicine. Part 4. Reference procedure for the measurement of catalytic concentration of alanine aminotransferase[J]. Clin Chem Lab Med, 2002, 40:718-724.

Symbol Interpretation

IVD	In Vitro Diagnostic Medical Device	LOT	Batch Code
	Consult Instructions for Use		Use-By Date
REF	Catalogue Number		Manufacturer
	Temperature Limit		Date of Manufacture
CE	CE marking of conformity	EC REP	Authorized Representative in the European Community



Zybio Inc.
Floor 1 to Floor 5, Building 30, No. 6 of Taikang Road, Block C of Jianqiao Industrial Park, Dadukou District, 400082 Chongqing, PEOPLE'S REPUBLIC OF CHINA.

Web: www.zybio.com

E-mail: info@zybio.com

Tel: +86 (0)23 6865 5509

Fax: +86 (0)23 6869 9779



Lotus NL B.V.
Koningin Julianaplein 10, 1e Verd, 2595AA, The Hague, Netherlands.

Current Version: 02 Date of Issue: May, 2022

Instructions for Use of High Density Lipoprotein Cholesterol (HDL-C) Kit (Enzymatic Method)

Package Specification

REF	Reagent	Systems
01.09.02.03.EC.01	R1 30 mL x 3 R2 10 mL x 3 Calibrator 1 Level x 1.0 mL x 1 Control 2 Levels x 1.0 mL x 1	Zybio EXC200/220
01.09.02.03.EC.02	R1 45 mL x 2 R2 15 mL x 2 Calibrator 1 Level x 1.0 mL x 1 Control 2 Levels x 1.0 mL x 1	Hitachi 7180 Zybio EXC400/420

Intended Use

In vitro test for the quantitative determination of high density lipoprotein cholesterol (HDL-C) concentration in human samples (serum). Clinically, it is mainly used for diagnosing hypercholesterolemia, coronary heart disease, and atherosclerosis.

Summary

A high HDL-C protects against CHD, as HDL-C is responsible for removing cholesterol from the periphery to the liver for catabolism. In the United States, the average woman has an HDL-C of 55 mg/dL and the average man has an HDL-C of 45 mg/dL. A 1-mg/dL increase in HDL-C has generally been associated with a 2% decrease in risk of CHD. The NCEP ATP III considers an HDL-C > 60 mg/dL as a negative risk factor, whereas HDL-C < 40 mg/dL is considered a positive risk factor for CHD. However, a study of postmenopausal women with CHD found that 20% had "protective" HDL-C of > 60 mg/dL. These women actually had fewer CHD risk factors than other women in the cohort. Because only 20% of women in the study had a high HDL-C vs 30% of women in the general population, this higher level of HDL-C did have some apparent role in CHD prevention. However, as an HDL-C of > 60 mg/dL is the 85th percentile for men vs the 70th percentile in women, it seems probable that the level of HDL-C considered cardioprotective may need to be higher in women. An HDL-C of 70 mg/dL is the 85th percentile for women, and this may be a more appropriate designated "negative risk factor" than an HDL-C of 60 mg/dL in women.

Principle

The kit uses enzymatic method to determine the concentration of high density lipoprotein cholesterol (HDL-C) in samples.

The non-HDL components of serum, such as chyle (CM), low-density lipoprotein (LDL) and very low-density lipoprotein (VLDL), are consumed by hydrolysis reactions under the action of enzymes such as cholesterol oxidase and peroxidase in reagent R1. High-density lipoprotein cholesterol was released under the action of denaturant in reagent R2, and HDL-C was detected by the following reactions.

- Cholesteryl Ester + H₂O $\xrightarrow{\text{Cholesterol Esterase}}$ Cholesterol + Fatty acid
- Cholesterol + O₂ $\xrightarrow{\text{Cholesterol Oxidase}}$ Cholestenone + H₂O₂
- 2H₂O₂ + 4-AAP + TOOS $\xrightarrow{\text{Peroxidase}}$ Quinonimine + 4H₂O

The production of quinonimine products causes an increase in absorbance at 546 nm, which is directly proportional to the HDL-C concentration in the sample. The concentration of HDL-C in the sample shall be calculated by measuring the absorbance change at 546 nm and comparing with that in calibrator treated in the same manner.

Reagents Components and Concentration

Components	Main Constituents	Concentration
R1	Cholesterol Oxidase	0.8-1.2 kU /L
	4-Aminoantipyrine (4-AAP)	0.05-0.15 g/L
	Peroxidase	1.6-2.4 kU /L

	2-Morpholinoethane sulfonic Acid (MES)	5.5-7.5 g/L
R2	Sodium 3-(N-ethyl-3-methylanilino)-2-hydroxypropanesulfonate (TOOS)	0.5-1.5 g/L
	Cholesterol Esterase	0.8-1.2 kU/L
	2-Morpholinoethane sulfonic Acid (MES)	5.5-7.5 g/L
Calibrator	Bovine Serum	Refer to the label for marked value of HDL-C concentration
	Sucrose	35 g/L
Control	Bovine Serum	Refer to the label for marked value of HDL-C concentration
	Sucrose	35 g/L

The components in different batches are non-interchangeable.
The measurement system can be traceable to JCCRM 224-16.
The target value of control has batch specificity.

Storage and Validity

- The reagents should be stored at 2 - 8 °C and kept away from freezing. The unopened reagents are valid for 18 months. Summer transportation with attention to refrigeration.
- Once opened, the reagents are stable for 1 month at 2 - 8 °C. For reagents not in use, the cap should be tightened to avoid contamination.
- To ensure accuracy, calibrator and control are stored at 2 - 8 °C after reconstitution and used only on the same day.
- The production date and expiration date are available on package insert.

System Information

Hitachi 7180, Zybio EXC400/420, Zybio EXC200/220 Chemistry Analyzer. Other models can be used after verification.

Specimen Information

Serum is suitable for samples, which shall be separated in time after collection to avoid hemolysis. Samples are stable for 6 days at 2 - 8 °C and 21 days at - 20 °C. Avoid repeated freezing and thawing.

Warnings and Precautions

- For in vitro diagnostic use only. In case of contact, immediately flush the affected area with plenty of water.
- The dosage of reagents and samples can be appropriately increased or decreased according to various instruments, under the condition that the volume proportion of reagents and samples is invariable.
- The necessary precautions should be taken to use the reagents. The bottle cap should be tightened immediately after use to avoid contamination.
- When reagent becomes turbid or the blank absorbance > 0.05, the reagent is failed and should be discarded.
- All samples and reaction wastes should be treated as sources of infection. The waste liquid generated during the experiment and the used packaging materials shall be collected and disposed in accordance with relevant regulations.
- Dedicated calibrator is recommended for use to ensure the accuracy of test values.
- The same sample tested with reagents from different manufacturers may lead to different measured values.
- Warning: This product contains ingredients of human or animal origin. At present, there is no way to completely guarantee that it is free from infectious substances and may be contaminated during use. This product and samples shall be considered as potential sources of infection, and operators should take protective measures and comply with laboratory safety regulations. All waste shall be disposed of in accordance with local regulations.

Test Process

1. Parameters

Method	End-Point Method	Sample/Reagent	3/320
Main Wavelength	546 nm	Reaction Temperature	37 °C
Sub Wavelength	700 nm	Reaction Time	10 min
Reaction Direction	+		

2. Operation

Addition	Blank	Calibration	Detection
Sample (μL)	/	/	3
Calibrator (μL)	/	3	/
Purified Water (μL)	3	/	/
Reagent 1 (μL)	240	240	240
Mix well, incubate at 37 °C for 5 min, and measure absorbance A ₁			
Reagent 2 (μL)	80	80	80
Mix well, incubate at 37 °C for 5 min, and measure absorbance A ₂ , calculate ΔA = A ₂ - A ₁ .			

3. Calibration

Use Zybio matched calibrator. Calibration cycle: Re-calibration is needed in case of replacement of the reagent batch, deviation of quality control, or the major repair and maintenance of equipment.

Calibrator reconstitution: Reconstitution with the amount of purified water labeled on the bottle accurately absorbed, leave for 30 minutes, and mix well before use.

4. Quality Control

To ensure the reliability of the test results, quality control is required after each calibration, and the quality control value should be within the specified range. If it is out of the specified range, the user should check the measurement system and recalibrate it if necessary. It is recommended that users conduct regular quality control according to their own needs to ensure the accuracy of the measurement system and the reliability of the test results.

Control reconstitution: Reconstitution with the amount of purified water labeled on the bottle accurately absorbed, leave for 30 minutes, and mix well before use.

5. Calculation

Linear calibration was used to draw the working curve. The concentration of high density lipoprotein cholesterol (HDL-C) in the sample can be calculated on the working curve based on its absorbance change value.

Reference Intervals

0.9~2.0 mmol/L

This reference interval is determined based on 95% distribution interval obtained from 200 healthy human specimens without related diseases, and is for reference only. It is recommended that each laboratory establish its own reference range.

Explanation of Results

If the concentration of HDL-C in the sample exceeds 4.00 mmol/L, the reduction mode of the chemistry analyzer will be used, or the high concentration sample shall be diluted with normal saline, with its report result multiplied by the dilution factor. The professional is responsible for the review of the test result, which may be affected by the subject's age, sex, or weight. The test values within the critical range should be confirmed by remeasuring, if it is obviously beyond the reference value range or if it is still beyond the reference range after confirmation detection, the target concentration in the sample is considered to be abnormal. If the test result is inconsistent with or even contrary to the clinical situation, cause analysis should be performed correspondingly.

Limitations

1. The deviation of test results caused by interferents is ≤ 10% if the concentrations of the following interferents are at or below the given values:

Substances	Concentrations
Hemoglobin	0.5 g/L
Chyle	0.30%
Bilirubin	300 μmol/L
Intralipid	1000 mg/dL

2. The test results of this product are only for clinical reference and cannot be used as the basis for diagnosis or exclusion of cases. The clinical diagnosis and treatment of patients should be combined with their symptoms/signs, medical history, other laboratory tests and treatment response. For achieve diagnostic purposes, the test results should be combined with clinical tests, medical history and other test results.

Performance Characteristics

- The reagent blank absorbance ≤ 0.05.
- Analytical sensitivity: at the test concentration of 1.00 mmol/L, the reagent absorbance change (ΔA) ≥ 0.04.
- Accuracy: relative deviation ≤ 10%.
- Precision: within-run CV ≤ 4%, between-run relative range ≤ 10%.
- Linear Range: [0.20, 4.00] mmol/L, the correlation coefficient (r) ≥ 0.995.
- [0.20, 1.00] mmol/L, the absolute deviation ≤ 0.10 mmol/L;
- [1.00, 4.00] mmol/L, the relative deviation ≤ 10%.
- Calibrator accuracy: relative deviation ≤ 10%.
- Calibrator homogeneity: between-vial CV ≤ 10%.
- Control accuracy: test value is within the allowable range of the marked value.
- Control homogeneity: between-vial CV ≤ 10%.

Materials Required (but not provided)

Chemistry analyzer, General lab equipment and consumable.

References

- [1] Gordon D, Probstfield J, Garrison R, et al. High-Density Lipoprotein Cholesterol and Cardiovascular Disease[J]. Am Heart Assoc, 1989, 79:8-15.
- [2] Warnick G, Nguyen R, Albers A. Comparison of improved precipitation methods for quantification of high density lipoprotein cholesterol[J]. Clin Chem, 1985, 31:217-222.
- [3] Warnick G, Cheung M, Albers J. Comparison of current methods for high-density lipoprotein cholesterol quantitation[J]. Clin Chem, 1979, 25:596-604.

Symbol Interpretation

IVD	In Vitro Diagnostic Medical Device	LOT	Batch Code
	Consult Instructions for Use		Use-By Date
REF	Catalogue Number		Manufacturer
	Temperature Limit		Date of Manufacture
CE	CE marking of conformity	EC REP	Authorized Representative in the European Community



Zybio Inc.
Floor 1 to Floor 5, Building 30, No. 6 of Taikang Road, Block C of Jianqiao Industrial Park, Dadukou District, 400082 Chongqing, PEOPLE'S REPUBLIC OF CHINA.
Web: www.zybio.com
E-mail: info@zybio.com
Tel: +86 (0)23 6865 5509
Fax: +86 (0)23 6869 9779



Lotus NL B.V.
Koningin Julianaplein 10, 1e Verd, 2595AA, The Hague, Netherlands.

Current Version: 02 Date of Issue: May, 2022

Instructions for Use of Lactate Dehydrogenase (LDH) Kit (Rate Method)

Package Specification

REF	Reagent	Systems
01.09.04.06.EC.01	R1 30 mL × 3	Zybio EXC200/220
	R2 7.5 mL × 3	
01.09.04.06.EC.02	R1 48 mL × 2	Hitachi 7180
	R2 12 mL × 2	Zybio EXC400/420

Intended Use


In vitro test for the quantitative determination of lactate dehydrogenase activity in human samples (serum). Clinically, it is mainly used as an aid to diagnosis of myocardial infarction and hepatopathy.

Summary

Lactate dehydrogenase is a kind of NAD-dependent kinase, which has three subunits, LDHA, LDHB and LDHC, and can constitute six tetrameric isoenzymes. Animal lactate dehydrogenase is a tetramer composed of 4 subunits, 5 LDH isozymes (LDH1-5) composed of common A and B subunits, and only one LDH isozyme (LDH-C4) composed of C subunit. Lactate dehydrogenase is a metalloprotein containing zinc ions, with a molecular weight of 135-140 kD. It is one of the important enzymes for anaerobic glycolysis and gluconeogenesis of sugars. It can catalyze the reduction and oxidation reaction between propionic acid and L-lactic acid, and can also catalyze the related α -keto acid. LDH is widely present in human tissues, with the highest content in the kidney, followed by the myocardium and bony muscle. LDH in red blood cells is about 100 times higher than in normal serum.

Lactate dehydrogenase is a key enzyme in microorganisms that catalyzes the production of benzolactic acid (also known as 2-hydroxy-3-phenylpropionic acid) from phenylpyruvate. Lactate dehydrogenase is a crucial oxidoreductase in the glycolytic pathway in organisms, which can reversibly catalyze the oxidation of lactate to pyruvate, and this catalytic reaction is the end product of anaerobic glycolysis. Lactate dehydrogenase is mainly found in animal tissues such as heart muscle, liver, kidney, skeletal muscle, or lung. Lactate dehydrogenase measurements are commonly used in the diagnosis of myocardial infarction, liver disease, and certain malignancies.

Principle

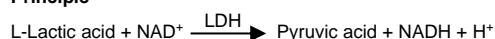

The activity of lactate dehydrogenase (LDH) in the sample can be detected by measuring the increase rate of the absorbance at 340 nm.

Reagents Components and Concentration

Components	Main Constituents	Concentration
R1	Lactate	>6 mmol/L
	Trometamol (Tris) buffer	100 mmol/L
R2	Nicotinamide adenine dinucleotide (NAD ⁺)	>12 mmol/L

The components in different batches are non-interchangeable.

Storage and Validity

- The reagents should be stored at 2 - 8 °C and kept away from direct light and freezing. The unopened reagents are valid for 12 months.
- Once opened, the reagents are stable for 30 days at 2 - 8 °C. For reagents not in

Floor 1 to Floor 5, Building 30, No. 6 of Taikang Road, Block C of Jianqiao Industrial Park, Dadukou District, 400082 Chongqing, PEOPLE'S REPUBLIC OF CHINA
Tel: +86 (0)23 6865 5509 Fax: +86 (0)23 6869 9779 E-mail: info@zybio.com Web: www.zybio.com

use, the cap should be tightened to avoid contamination.

- The production date and expiration date are available on package insert.

System Information

Hitachi 7180, Zybio EXC400/420, Zybio EXC200/220 Chemistry Analyzer. Other models can be used after verification.

Specimen Information

Serum is suitable for samples, which are stable for 3 days at 2 - 8 °C and for 30 days at - 20 °C. Avoid hemolysis and repeated freezing and thawing.

Warnings and Precautions

- For in vitro diagnostic use only. In case of contact, immediately flush the affected area with plenty of water.
- The dosage of reagents and samples can be appropriately increased or decreased according to various instruments, under the condition that the volume proportion of reagents and samples is invariable.
- The necessary precautions should be taken to use the reagents. The bottle cap should be tightened immediately after use to avoid contamination.
- When reagent becomes turbid or the blank absorbance > 0.500, the reagent is failed and should be discarded.
- All samples and reaction wastes should be treated as sources of infection. The waste liquid generated during the experiment and the used packaging materials shall be collected and disposed in accordance with relevant regulations.
- The same sample tested with reagents from different manufacturers may lead to different measured values.
- Warning: This product contains ingredients of human or animal origin. At present, there is no way to completely guarantee that it is free from infectious substances and may be contaminated during use. This product and samples shall be considered as potential sources of infection, and operators should take protective measures and comply with laboratory safety regulations. All waste shall be disposed of in accordance with local regulations.

Test Process

1. Parameters

Method	Rate Method	Sample/Reagent	1/50
Main Wavelength	340 nm	Reaction Temperature	37 °C
Sub Wavelength	405 nm	Reaction Time	10 min
Reaction Direction	+		

2. Operation

Addition	Blank	Calibration	Detection
Sample (μL)	/	/	5
Calibrator (μL)	/	5	/
Purified Water (μL)	5	/	/
Reagent 1 (μL)	200	200	200
Mix well, incubate at 37 °C for 5 min			
Reagent 2 (μL)	50	50	50
Mix well, after 2 min, measure the average absorbance change rate $\Delta A/\text{min}$ within 3 min.			

3. Calibration

Use Zybio Clinical Chemistry Multi-analyte Calibrator or Randox multi-analyte calibrator. Calibration cycle: Re-calibration is needed in case of replacement of the reagent batch, deviation of quality control, or the major repair and maintenance of equipment.

4. Quality Control

To ensure the reliability of the test results, quality control is required after each calibration, and the quality control value should be within the specified range. If it is out of the specified range, the user should check the measurement system and recalibrate it if necessary. It is recommended that users conduct regular quality control according to their own needs to ensure the accuracy of the measurement system and the reliability of the test results.

5. Calculation

Linear calibration was used to draw the working curve. The concentration of lactate dehydrogenase (LDH) in the sample can be calculated on the working curve based on its absorbance change rate.

Reference Intervals

105~245 U/L

This reference interval is determined based on 95% distribution interval obtained from 200 healthy human specimens without related diseases, and is for reference only. It is recommended that each laboratory establish its own reference range.

Explanation of Results

If the concentration of LDH in the sample exceeds 800 U/L, the reduction mode of the chemistry analyzer will be used, or the high concentration sample shall be diluted with normal saline, with its report result multiplied by the dilution factor.

The professional is responsible for the review of the test result, which may be affected by the subject's age, sex, or weight. The test values within the critical range should be confirmed by remeasuring, if it is obviously beyond the reference value range or if it is still beyond the reference range after confirmation detection, the target concentration in the sample is considered to be abnormal. If the test result is inconsistent with or even contrary to the clinical situation, cause analysis should be performed correspondingly.

Limitations

1. The deviation of test results caused by interferents is less than 10% if the concentrations of the following interferents are at or below the given values:

Substances	Concentrations
Vc	0.5 g/L
Chyle	0.30%
Bilirubin	342 μ mol/L

2. The test results of this product are only for clinical reference and cannot be used as the basis for diagnosis or exclusion of cases. The clinical diagnosis and treatment of patients should be combined with their symptoms/signs, medical history, other laboratory tests and treatment response. For achieve diagnostic purposes, the test results should be combined with clinical tests, medical history and other test results.

Performance Characteristics

1. The reagent blank absorbance ≤ 0.500 ; the reagent blank absorbance change rate ($\Delta A/\text{min}$) $\leq 0.002 A/\text{min}$.

2. Analytical sensitivity: at the test concentration of 200 U/L, the reagent absorbance change rate ($\Delta A/\text{min}$) ≥ 0.005 .

3. Accuracy: relative deviation $\leq 10\%$.

4. Precision: within-run CV $\leq 5\%$, between-run relative range $\leq 10\%$.

5. Linear Range:

[25, 800] U/L, the correlation coefficient (r) ≥ 0.990 .

[25, 100] U/L, the absolute deviation ≤ 10 U/L;

(100, 800] U/L, the relative deviation $\leq 10\%$.

Floor 1 to Floor 5, Building 30, No. 6 of Taikang Road, Block C of Jianqiao Industrial Park, Dadukou District, 400082 Chongqing, PEOPLE'S REPUBLIC OF CHINA
Tel: +86 (0)23 6865 5509 Fax: +86 (0)23 6869 9779 E-mail: info@zybio.com Web: www.zybio.com

Materials Required (but not provided)

Chemistry analyzer, Zybio Clinical Chemistry Multi-analyte Calibrator or Randox multi-analyte calibrator, Control, General lab equipment and consumable.

References

[1] Eigentler T, Figl A, Krex D, et al. Number of metastases, serum lactate dehydrogenase level, and type of treatment are prognostic factors in patients with brain metastases of malignant melanoma[J]. Cancer, 2011, 117:1697-1703.

Symbol Interpretation

IVD	In Vitro Diagnostic Medical Device	LOT	Batch Code
	Consult Instructions for Use		Use-By Date
REF	Catalogue Number		Manufacturer
	Temperature Limit		Date of Manufacture
CE	CE marking of conformity	EC REP	Authorized Representative in the European Community



Zybio Inc.

Floor 1 to Floor 5, Building 30, No. 6 of Taikang Road, Block C of Jianqiao Industrial Park, Dadukou District, 400082 Chongqing, PEOPLE'S REPUBLIC OF CHINA.

Web: www.zybio.com

E-mail: info@zybio.com

Tel: +86 (0)23 6865 5509

Fax: +86 (0)23 6869 9779



Lotus NL B.V.

Koningin Julianaplein 10, 1e Verd, 2595AA, The Hague, Netherlands.

Current Version: 02 Date of Issue: May, 2022

Instructions for Use of Low Density Lipoprotein Cholesterol (LDL-C) Kit (Enzymatic Method)

Package Specification

REF	Reagent	Systems
01.09.02.00.EC.01	R1 30 mL × 3 R2 10 mL × 3 Calibrator 1 Level × 1.0 mL × 1 Control 2 Levels × 1.0 mL × 1	Zybio EXC200/220
01.09.02.00.EC.02	R1 45 mL × 2 R2 15 mL × 2 Calibrator 1 Level × 1.0 mL × 1 Control 2 Levels × 1.0 mL × 1	Hitachi 7180 Zybio EXC400/420

Intended Use


In vitro test for the quantitative determination of low density lipoprotein cholesterol (LDL-C) concentration in human samples (serum). Clinically, it is mainly used as an aid to diagnosis of hypercholesterolemia, coronary heart disease, and atherosclerosis.

Summary


Low Density Lipoprotein (LDL) play a key role in causing and influencing the progression of atherosclerosis and, in particular, coronary sclerosis. The LDLs are derived from VLDLs (Very Low Density Lipoproteins) rich in triglycerides by the action of various lipolytic enzymes and are synthesized in the liver. The elimination of LDL from plasma takes place mainly by liver parenchymal cells via specific LDL receptors. Elevated LDL concentrations in blood and an increase in their residence time coupled with an increase in the biological modification rate results in the destruction of the endothelial function and a higher LDL-cholesterol uptake in the monocyte/macrophage system as well as by smooth muscle cells in vessel walls. The majority of cholesterol stored in atherosclerotic plaques originates from LDL. The LDL-cholesterol value is the most powerful clinical predictor among all of the single parameters with respect to coronary atherosclerosis. Therefore, therapies focusing on lipid reduction primarily target the reduction of LDL-cholesterol which is then expressed in an improvement of the endothelial function, prevention of atherosclerosis and reducing its progression as well as preventing plaque rupture.

Principle

1. The surface-active ingredients in R1 inhibit low density lipoprotein in serum, while high-density lipoprotein and very low-density lipoprotein are consumed by the reaction catalyzed by cholesterol enzyme.

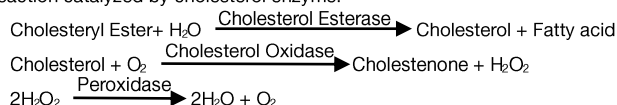

2. LDL-C is only measured after the surfactant in R2 has released low-density lipoprotein.

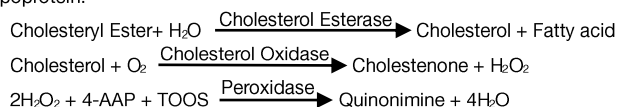

3. The absorbance of quinonimine is directly proportional to the content of cholesterol. The content of low-density lipoprotein cholesterol (LDL-C) in the sample can be calculated by measuring the absorbance change value at 546 nm.

Reagents Components and Concentration

Components	Main Constituents	Concentration
R1	Ascorbate Oxidase	> 3000 U/L
	Cholesterol Oxidase	> 400 U/L
	Cholesterol Esterase	> 500 U/L
	Sodium 3-(N-ethyl-3-methylanilino)-2-hydroxypropanesulfonate (TOOS)	0.8 mmol/L

	1,4-Piperazinebis (ethanesulfonic acid) buffer (PIPES)	100 mmol/L
R2	4-Aminoantipyrine (4-AAP)	4 mmol/L
	1,4-Piperazinebis (ethanesulfonic acid) buffer (PIPES)	100 mmol/L
	Peroxidase	> 1500 U/L
	Surfactant	Appropriate amount
Calibrator	Sucrose	10 g/L
	Bovine Serum	Refer to the label for marked value of LDL-C concentration
Control	Sucrose	10 g/L
	Bovine Serum	Refer to the label for marked value of LDL-C concentration

The measurement system can be traceable to JCCRM 224-16.

The components in different batches are non-interchangeable.

The target value of control has batch specificity.

Storage and Validity

1. The reagents should be stored at 2 - 8 °C and kept away from direct light and freezing. The unopened reagents are valid for 18 months. Summer transportation with attention to refrigeration.
2. Once opened, the reagents are stable for 1 month at 2 - 8 °C. For reagents not in use, the cap should be tightened to avoid contamination.
3. To ensure accuracy, calibrator and control are stored at 2 - 8 °C after reconstitution and used only on the same day.
4. The production date and expiration date are available on package insert.

System Information

Hitachi 7180, Zybio EXC400/420, Zybio EXC200/220 Chemistry Analyzer. Other models can be used after verification.

Specimen Information

Serum is suitable for samples, which shall be separated in time after collection to avoid hemolysis. Samples are stable for 6 days at 2 - 8 °C and 3 weeks at - 20 °C. Avoid repeated freezing and thawing.

Warnings and Precautions

1. For in vitro diagnostic use only. In case of contact, immediately flush the affected area with plenty of water.
2. The dosage of reagents and samples can be appropriately increased or decreased according to various instruments, under the condition that the volume proportion of reagents and samples is invariable.
3. The necessary precautions should be taken to use the reagents. The bottle cap should be tightened immediately after use to avoid contamination.
4. When reagent becomes turbid or the blank absorbance > 0.05, the reagent is failed and should be discarded.
5. All samples and reaction wastes should be treated as sources of infection. The waste liquid generated during the experiment and the used packaging materials shall be collected and disposed in accordance with relevant regulations.
6. Dedicated calibrator is recommended for use to ensure the accuracy of test values.
7. The same sample tested with reagents from different manufacturers may lead to different measured values.
8. Warning: This product contains ingredients of human or animal origin. At present, there is no way to completely guarantee that it is free from infectious substances and may be contaminated during use. This product and samples shall be considered as potential sources of infection, and operators should take protective measures and comply with laboratory safety regulations. All waste shall be disposed of in accordance with local regulations.

Test Process

1. Parameters

Method	End-Point Method	Sample/Reagent	3/320
Main Wavelength	546 nm	Reaction Temperature	37 °C
Sub Wavelength	700 nm	Reaction Time	10 min
Reaction Direction	+		

2. Operation

Addition	Blank	Calibration	Detection
Sample (μL)	/	/	3
Calibrator (μL)	/	3	/
Purified Water (μL)	3	/	/
Reagent 1 (μL)	240	240	240
Mix well, incubate at 37 °C for 5 min, and measure absorbance A_1			
Reagent 2 (μL)	80	80	80
Mix well, measure absorbance A_2 after 5 min, calculate $\Delta A = A_2 - A_1$			

3. Calibration

Use Zybio matched calibrator. Calibration cycle: Re-calibration is needed in case of replacement of the reagent batch, deviation of quality control, or the major repair and maintenance of equipment.

Calibrator reconstitution: Reconstitution with the amount of purified water labeled on the bottle accurately absorbed, leave for 30 minutes, and mix well before use.

4. Quality Control

To ensure the reliability of the test results, quality control is required after each calibration, and the quality control value should be within the specified range. If it is out of the specified range, the user should check the measurement system and recalibrate it if necessary. It is recommended that users conduct regular quality control according to their own needs to ensure the accuracy of the measurement system and the reliability of the test results.

Control reconstitution: Reconstitution with the amount of purified water labeled on the bottle accurately absorbed, leave for 30 minutes, and mix well before use.

5. Calculation

Linear calibration was used to draw the working curve. The concentration of low density lipoprotein cholesterol (LDL-C) in the sample can be calculated on the working curve based on its absorbance change value.

Reference Intervals

≤3.36 mmol/L

This reference interval is determined based on 95% distribution interval obtained from 200 healthy human specimens without related diseases, and is for reference only. It is recommended that each laboratory establish its own reference range.

Explanation of Results

If the concentration of LDL-C in the sample exceeds 11.60 mmol/L, the reduction mode of the chemistry analyzer will be used, or the high concentration sample shall be diluted with normal saline, with its report result multiplied by the dilution factor. The professional is responsible for the review of the test result, which may be affected by the subject's age, sex, or weight. The test values within the critical range should be confirmed by remeasuring, if it is obviously beyond the reference value range or if it is still beyond the reference range after confirmation detection, the target concentration in the sample is considered to be abnormal. If the test result is inconsistent with or even contrary to the clinical situation, cause analysis should be performed correspondingly.

Limitations

1. The deviation of test results caused by interferents is ≤ 10% if the concentrations of the following interferents are at or below the given values:

Substances	Concentrations
Vc	0.5 g/L
Hemoglobin	5 g/L
Chyle	0.30%
Bilirubin	342 μmol/L
Intralipid	1000 mg/dL

2. The test results of this product are only for clinical reference and cannot be used as the basis for diagnosis or exclusion of cases. The clinical diagnosis and treatment of patients should be combined with their symptoms/signs, medical history, other laboratory tests and treatment response. For achieve diagnostic purposes, the test results should be combined with clinical tests, medical history and other test results.

Performance Characteristics

- The reagent blank absorbance ≤ 0.05.
- Analytical sensitivity: at the test concentration of 1.00 mmol/L, the reagent absorbance change (ΔA) ≥ 0.03.
- Accuracy: relative deviation ≤ 10%.
- Precision: within-run CV ≤ 3%, between-run relative range ≤ 10%.
- Linear Range: [0.20, 11.60] mmol/L, the correlation coefficient (r) ≥ 0.995. [0.20, 3.00] mmol/L, the absolute deviation ≤ 0.30 mmol/L; [3.00, 11.60] mmol/L, the relative deviation ≤ 10%.
- Calibrator accuracy: relative deviation ≤ 10%.
- Calibrator homogeneity: between-vial CV ≤ 10%.
- Control accuracy: test value is within the allowable range of the marked value.
- Control homogeneity: between-vial CV ≤ 10%.

Materials Required (but not provided)

Chemistry analyzer, General lab equipment and consumable.

References

[1] Davidson M. Low-density lipoprotein cholesterol, non-high-density lipoprotein, apolipoprotein, or low-density lipoprotein particle: what should clinicians measure? [J]. J Am Coll Cardiol, 2012, 60:2616-2617.

Symbol Interpretation

IVD	In Vitro Diagnostic Medical Device	LOT	Batch Code
	Consult Instructions for Use		Use-By Date
REF	Catalogue Number		Manufacturer
	Temperature Limit		Date of Manufacture
CE	CE marking of conformity	EC REP	Authorized Representative in the European Community



Zybio Inc.
Floor 1 to Floor 5, Building 30, No. 6 of Taikang Road, Block C of Jianqiao Industrial Park, Dadukou District, 400082 Chongqing, PEOPLE'S REPUBLIC OF CHINA.
Web: www.zybio.com
E-mail: info@zybio.com
Tel: +86 (0)23 6865 5509
Fax: +86 (0)23 6869 9779



Lotus NL B.V.
Koningin Julianaplein 10, 1e Verd, 2595AA, The Hague, Netherlands.

Current Version: 02 Date of Issue: May, 2022

Instructions for Use of Urea (UREA) Kit (Urease-GLDH Method)

Package Specification

REF	Reagent	Systems
01.09.01.06.EC.01	R1 30 mL x 3	Zybio EXC200/220
	R2 7.5 mL x 3	
01.09.01.06.EC.02	R1 48 mL x 2	Hitachi 7180
	R2 12 mL x 2	Zybio EXC400/420



Intended Use

In vitro test for the quantitative determination of urea concentration in human samples (serum or plasma). Clinically, it is mainly used as one of the evaluation indicators of renal function.

Summary

Urea is the major end product of protein nitrogen metabolism. It is synthesized by the urea cycle in the liver from ammonia which is produced by amino acid deamination. Urea is excreted mostly by the kidneys but minimal amounts are also excreted in sweat and degraded in the intestines by bacterial action. Determination of blood urea nitrogen is the most widely used screening test for renal function. When used in conjunction with serum creatinine determinations it can aid in the differential diagnosis of the three types of azotemia: prerenal, renal and postrenal. Elevations in blood urea nitrogen concentration are seen in inadequate renal perfusion, shock, diminished blood volume (prerenal causes), chronic nephritis, nephrosclerosis, tubular necrosis, glomerular nephritis (renal causes) and urinary tract obstruction (postrenal causes). Transient elevations may also be seen during periods of high protein intake. Unpredictable levels occur with liver diseases.

Principle

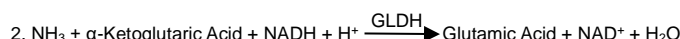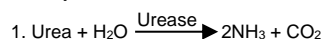

Oxidation of NADH to NAD⁺ causes a decrease in absorbance at 340 nm, which is directly proportional to the Urea concentration in the sample.

Reagents Components and Concentration

Components	Main Constituents	Concentration
R1	Trometamol (Tris) buffer	100 mmol/L
	Nicotinamide adenine dinucleotide (NADH)	0.3 mmol/L
R2	α-Ketoglutaric Acid	10 mmol/L
	Urease	6.0 kU/L
	Glutamate dehydrogenase (GLDH)	2.0 kU/L

The components in different batches are non-interchangeable.

Storage and Validity

- The reagents should be stored at 2 - 8 °C and kept away from direct light and freezing. The unopened reagents are valid for 12 months.
- Once opened, the reagents are stable for 30 days at 2 - 8 °C. For reagents not in use, the cap should be tightened to avoid contamination.
- The production date and expiration date are available on package insert.

System Information

Hitachi 7180, Zybio EXC400/420, Zybio EXC200/220 Chemistry Analyzer. Other models can be used after verification.

Specimen Information

Serum or plasma (heparin or EDTA anticoagulation) is suitable for samples, which are stable for 3 days at 2 - 8 °C and for 30 days at - 20 °C. Avoid repeated freezing and thawing.

Warnings and Precautions

- For in vitro diagnostic use only. In case of contact, immediately flush the affected area with plenty of water.
- The dosage of reagents and samples can be appropriately increased or decreased according to various instruments, under the condition that the volume proportion of reagents and samples is invariable.
- The necessary precautions should be taken to use the reagents. The bottle cap should be tightened immediately after use to avoid contamination.
- When reagent becomes turbid or the blank absorbance < 1.000, the reagent is failed and should be discarded.
- All samples and reaction wastes should be treated as sources of infection. The waste liquid generated during the experiment and the used packaging materials shall be collected and disposed in accordance with relevant regulations.
- The same sample tested with reagents from different manufacturers may lead to different measured values.
- Warning: This product contains ingredients of human or animal origin. At present, there is no way to completely guarantee that it is free from infectious substances and may be contaminated during use. This product and samples shall be considered as potential sources of infection, and operators should take protective measures and comply with laboratory safety regulations. All waste shall be disposed of in accordance with local regulations.

Test Process

1. Parameters

Method	Rate Method	Sample/Reagent	1/100
Main Wavelength	340 nm	Reaction Temperature	37 °C
Sub Wavelength	405 nm	Reaction Time	10 min
Reaction Direction	-		

2. Operation

Addition	Blank	Calibration	Detection
Sample (μL)	/	/	3
Calibrator (μL)	/	3	/
Purified Water (μL)	3	/	/
Reagent 1 (μL)	240	240	240
Mix well, incubate at 37 °C for 5 min			
Reagent 2 (μL)	60	60	60
Mix well, after 1 min, measure the absorbance change within 2 min, and calculate the absorbance change rate ΔA/ min.			

3. Calibration

Use Zybio Clinical Chemistry Multi-analyte Calibrator or Randox multi-analyte calibrator. Calibration cycle: Re-calibration is needed in case of replacement of the reagent batch, deviation of quality control, or the major repair and maintenance of equipment.

4. Quality Control

To ensure the reliability of the test results, quality control is required after each calibration, and the quality control value should be within the specified range. If it is out of the specified range, the user should check the measurement system and recalibrate it if necessary. It is recommended that users conduct regular quality control according to their own needs to ensure the accuracy of the measurement system and the reliability of the test results.

5. Calculation

Linear calibration was used to draw the working curve. The concentration of urea (UREA) in the sample can be calculated on the working curve based on its absorbance change rate.

Reference Intervals

1.7~8.3 mmol/L (10~50 mg/dL)

This reference interval is determined based on 95% distribution interval obtained from 200 healthy human specimens without related diseases, and is for reference only. It is recommended that each laboratory establish its own reference range.

Explanation of Results

If the concentration of UREA in the sample exceeds 40.0 mmol/L, the reduction mode of the chemistry analyzer will be used, or the high concentration sample shall be diluted with normal saline, with its report result multiplied by the dilution factor.

The professional is responsible for the review of the test result, which may be affected by the subject's age, sex, or weight. The test values within the critical range should be confirmed by remeasuring, if it is obviously beyond the reference value range or if it is still beyond the reference range after confirmation detection, the target concentration in the sample is considered to be abnormal. If the test result is inconsistent with or even contrary to the clinical situation, cause analysis should be performed correspondingly.

Limitations

1. The deviation of test results caused by interferents is less than 10% if the concentrations of the following interferents are at or below the given values:

Substances	Concentrations
Vc	0.5 g/L
Hemoglobin	5 g/L
Chyle	0.30%
Bilirubin	342 μmol/L

2. The test results of this product are only for clinical reference and cannot be used as the basis for diagnosis or exclusion of cases. The clinical diagnosis and treatment of patients should be combined with their symptoms/signs, medical history, other laboratory tests and treatment response. For achieve diagnostic purposes, the test results should be combined with clinical tests, medical history and other test results.

Performance Characteristics

1. The reagent blank absorbance ≥ 1.000 ; the reagent blank absorbance change rate ($\Delta A/\text{min}$) ≤ 0.04 .

2. Analytical sensitivity: at the test concentration of 7.5 mmol/L, the reagent absorbance change rate ($\Delta A/\text{min}$) ≥ 0.008 .

3. Accuracy: relative deviation $\leq 10\%$.

4. Precision: within-run CV $\leq 5\%$, between-run relative range $\leq 6\%$.

5. Linear Range:

[0.5, 40.0] mmol/L, the correlation coefficient (r) ≥ 0.990 .

[0.5, 5.0] mmol/L, the absolute deviation ≤ 0.5 mmol/L;

(5.0, 40.0] mmol/L, the relative deviation $\leq 10\%$.

Floor 1 to Floor 5, Building 30, No. 6 of Taikang Road, Block C of Jianqiao Industrial Park, Dadukou District, 400082 Chongqing, PEOPLE'S REPUBLIC OF CHINA
Tel: +86 (0)23 6865 5509 Fax: +86 (0)23 6869 9779 E-mail: info@zybio.com Web: www.zybio.com

Materials Required (but not provided)

Chemistry analyzer, Zybio Clinical Chemistry Multi-analyte Calibrator or Randox multi-analyte calibrator, Control, General lab equipment and consumable.

References

[1] Ai H, Chen K. Diagnostic Value of Blood Urea Nitrogen and Serum Creatinine in the Diagnosis of Early Diabetic Nephropathy[J]. Journal of Practical Medical Techniques, 2008, 15:431-433.

Symbol Interpretation

IVD	In Vitro Diagnostic Medical Device	LOT	Batch Code
	Consult Instructions for Use		Use-By Date
REF	Catalogue Number		Manufacturer
	Temperature Limit		Date of Manufacture
CE	CE marking of conformity	EC REP	Authorized Representative in the European Community



Zybio Inc.

Floor 1 to Floor 5, Building 30, No. 6 of Taikang Road, Block C of Jianqiao Industrial Park, Dadukou District, 400082 Chongqing, PEOPLE'S REPUBLIC OF CHINA.

Web: www.zybio.com

E-mail: info@zybio.com

Tel: +86 (0)23 6865 5509

Fax: +86 (0)23 6869 9779



Lotus NL B.V.

Koningin Julianaplein 10, 1e Verd, 2595AA, The Hague, Netherlands.

Current Version: 02 Date of Issue: May, 2022