# Mission® Urinalysis Reagent Strips and Urine Analyzers

Obtain reliable and cost-effective results with *Mission®* Urinalysis Reagent Strips and Urine Analyzers!

- *Accurate*
- *Reliable*
- *Convenient*

**ACON®**

Global Diagnostics for Local Markets™

# U120 Urine Analyzer

## Accurate

- Up to 120 tests/hour in Continuous Test Option
- Capable of reading 1 strip at a time in Single Test Option
- Test modes include Routine, STAT and QC
- Automatic calibration for accurate results and easy operation

## Reliable

- Can read up to 4 Strip combinations with 8, 9, 10, 11 parameters, additional strips with 1-11 parameters available upon request
- Minimal training required

## Convenient Operation

- Saves and recalls the last 2,000 results automatically
- Audible beep signals operator to dip strips in urine
- Can print up to 3 copies per test for convenient reviewing and easy record keeping
- Option to print results on sticker paper for quick and simple record management

## Easy Data Management

- Includes RS232C port for easy data transfer to an external computer or LIS
- Optional Barcode Reader to record patient ID

## Unique Lockout Functions *new!*

- Strip Lockout
  - Prevents using strips of another brand on the U120 Urine Analyzer
  - Requires barcode reader scan or manual entry of the canister code
- User Lockout
  - Eliminates unapproved users from testing
  - Up to 10 lab operators can perform testing, but only the lab administrator can change analyzer settings
- QC Lockout
  - Prevents testing without passing QC
  - QC tests can be performed once every 8 hours, day, week or month
  - Analyzer will alert when to run QC test
  - If QC tests fail, analyzer will switch to STAT mode and list "E" at the end of each test number

## Specifications

Feature	Specifications
Analyzer Type	Manual
Methodology	Reflectance Photometry
Detection	Photosensitive Diode
Throughput	Single Test Option: 60 tests/hour Continuous Test Option: 120 tests/hour
Test Modes	Routine, STAT and QC
Lockout Functions	Strip Lockout: Available Upon Request; User/QC Lockout: Included with option to turn ON/OFF
Memory	Last 2,000 results
Strip Incubation Time	1 Minute
Wavelength of Monochromatic LED	525 nm and 635 nm
Standard Strips	8, 9, 10, 11 Parameters (5 mm x 108 mm)
Additional Strips Available	1-11 Parameters (5 mm x 108 mm); see URS Parameters
Total Combinations Per Analyzer	4 Combinations
Analyzer Ports	Standard RS232C Port for Barcode Reader or Data Transfer USB Port for Data Transfer 25 Pin Parallel Port for External Printer
Capabilities	Internal Thermal Printer (included)      RS232C Barcode Reader (optional) Optional External Printer (not included)      USB or RS232C Data Transfer Cable (optional)
Major Readable Barcodes	Code 128, Code 39, Codabar (NW-7), Interleaved 25, UPC-A, UPC-E, EAN 8, EAN 13
Calibration	Automatic
Available Languages on the Screen	English and additional language(s)
Operating Conditions	0-40°C (32-104°F); ≤ 85% RH
Storage Conditions	-5-50°C (23-122°F); ≤ 90% RH
Power Source	100-240 VAC, 50-60 Hz
Dimensions (L x W x H)	27.2 cm x 26.9 cm x 14.6 cm (10.7" x 10.6" x 5.7")
Display Dimensions (L x W)	10.8 cm x 5.7 cm (4.2" x 2.2")
Weight	2.6 kg (5.7 lbs)

## Ordering Information

Product Name	Catalog No.	Components	Kit Box Dimensions (L x W x H) & Weight	Carton Dimensions (L x W x H) & Weight	Number of Kits/Carton
U120 Urine Analyzer	U111-101 <sup>†</sup>	1 Urine Analyzer 1 Strip holder 2 Printer Paper Rolls	2 Fuses (2.0A) 1 Power Cord 1 Quick Start Guide 1 Instruction Manual	42.0 cm x 41.5 cm x 31 cm; 5.0 kg 16.4" x 16.2" x 12.1"; 176.4 oz	1
U120 Urine Analyzer with Barcode Reader	U111-111 <sup>†</sup>	1 Urine Analyzer 1 Strip holder 2 Printer Paper Rolls 1 Barcode Reader (RS232C)	2 Fuses (2.0A) 1 Power Cord 1 Serial Splitter Cable (RS232C) 1 Quick Start Guide 1 Instruction Manual	44.5 cm x 44.5 cm x 40.0 cm; 5.5 kg 17.5" x 17.5" x 15.7"; 194 oz	1
Barcode Reader	U221-111 <sup>†</sup>	1 Barcode Reader (RS232C)	1 Serial Splitter Cable (RS232C)	23.6 cm x 10.8 cm x 7.8 cm; 0.482 kg 9.3" x 4.3" x 3.1"; 17.0 oz 24.8" x 14.6" x 11.8"; 423.3 oz	22
Printer Paper Rolls	U121-101	4 Printer Paper Rolls	Thermal Paper (0.06 m x 20 m): 200 results/roll Sticker Paper (0.06 m x 9 m): 100 results/roll	12.0 cm x 12.0 cm x 6.5 cm; 0.36 kg 4.7" x 4.7" x 2.6"; 12.7 oz 63.0 cm x 37.0 cm x 30.0 cm; 19.4 kg 24.8" x 14.6" x 11.8"; 684.3 oz 12.0 cm x 12.0 cm x 6.5 cm; 0.4 kg 4.7" x 4.7" x 2.6"; 14.1 oz 63.0 cm x 37.0 cm x 30.0 cm; 21.4 kg 24.8" x 14.6" x 11.8"; 684.3 oz; 754.9 oz	50
U120 Data Transfer Kit	U221-131 <sup>†</sup>	1 Data Transfer Cable (RS232C)	1 Package Insert	16.0 cm x 13.0 cm x 3.5 cm; 0.147 kg 6.3" x 5.1" x 1.4"; 5.2 oz 25.0 cm x 21.0 cm x 15.0 cm; 1.36 kg 9.8" x 8.3" x 5.9"; 48.0 oz	8

✓ CE Marked for sale in the European Community **CE**  
 † Cleared for US 510(k)

# U500 Urine Analyzer

## Accurate and Efficient

- Up to 500 tests/hour for medium/large volume sample testing
- Professional accuracy equivalent to market leader
- Automatic strip detection and alignment for better efficiency
- Test modes include Routine, STAT and QC

## Easy to Operate

- Large touch screen LCD offers simple menu navigation
- Uniquely designed strip platform/waste tray unit for easy one-step cleaning

## Convenient

- Automatic calibration and waste disposal reduce hands-on time
- Can read strips with 8, 9, 10, 11 parameters, additional strips with 1-11 parameters available upon request
- Strip selection of up to 4 combinations for analyzer reading
- Stores up to 2,000 records and automatically flags abnormal results
- Capable of printing results on sticker paper for quick and easy record management

## Data Management Capability

- Includes RS232C port for easy data transfer to an external computer or LIS
- Optional Barcode Reader to record patient ID

## Unique Lockout Functions Coming Soon!

- Strip Lockout
  - Prevents using strips of another brand on the U500 Urine Analyzer
  - Requires barcode reader scan or manual entry of the canister code
- User Lockout
  - Eliminates unapproved users from testing
  - Up to 10 lab operators can perform testing, but only the lab administrator can change analyzer settings
- QC Lockout
  - Prevents testing without passing QC
  - QC tests can be performed once every 8 hours, day, week or month
  - Analyzer will alert when to run QC test
  - If QC tests fail, analyzer will switch to STAT mode and list "E" at the end of each test number

## Specifications

Feature	Specifications
Analyzer Type	Semi-Automatic
Methodology	Reflectance Photometry
Detection	Photosensitive Diode
Throughput	500 tests/hour (Measuring cycle: 7 seconds/test)
Test Modes	Routine, STAT and QC
Lockout Functions	Strip Lockout: Available Upon Request; User/QC Lockout: Included with option to turn ON/OFF
Memory	Last 2,000 Records
Strip Incubation Time	1 Minute
Wavelength	525 and 635 nm
Standard Strips	8, 9, 10, 11 Parameters (5 mm x 108 mm)
Additional Strips Available	1-11 Parameters (5 mm x 108 mm); see URS Parameters
Total Combinations Per Analyzer	4 Combinations
Waste Disposal Capacity	Up to 150 Strips
Analyzer Ports	Standard RS232C Port for Barcode Reader or Data Transfer 25 Pin Parallel Port for External Printer
Capabilities	Internal Thermal Printer (included)      RS232C Barcode Reader (optional) Optional External Printer (not included)      RS232C Data Transfer Cable (optional)
Major Readable Barcodes	Code 128, Code 39, Codabar (NW-7), Interleaved 25, UPC-A, UPC-E, EAN 8, EAN 13
Calibration	Automatic
Available Languages on the Screen	English and additional language(s)
Operating Conditions	0-40°C (32-104°F); ≤85% RH
Storage Conditions	-5-50°C (23-122°F); ≤90% RH
Power Source	100-240 VAC, 50-60 Hz
Dimensions (L x W x H)	36.6 cm x 28.3 cm x 19.5cm (14.4" x 11.1" x 7.7")
Display Dimensions (L x W)	11.5 cm x 9.0 cm (4.5" x 3.5")
Weight	4.0 kg (8.8 lbs)

## Ordering Information

Product Name	Catalog No.	Components	Kit Box Dimensions (L x W x H) & Weight	Carton Dimensions (L x W x H) & Weight	Number of Kits/Carton
U500 Urine Analyzer	U211-101 <sup>✓</sup>	1 Urine Analyzer 1 Strip Platform/Waste Tray 2 Printer Paper Rolls	2 Fuses (2.0A) 1 Power Cord 1 Instruction Manual	51.0 cm x 42.0 cm x 38.5 cm; 7 kg 20.1" x 16.5" x 15.2"; 246.9 oz	1
U500 Urine Analyzer with Barcode Reader	U211-111 <sup>✓</sup>	1 Urine Analyzer 1 Strip Platform/Waste Tray 2 Printer Paper Rolls 1 Barcode Reader (RS232C)	2 Fuses (2.0A) 1 Power Cord 1 Serial Splitter Cable (RS232C) 1 Instruction Manual	55.0 cm x 55.0 cm x 55.0cm; 9.2 kg 21.7" x 21.7" x 21.7"; 324.5 oz	1
Barcode Reader	U221-111 <sup>✓†</sup>	1 Barcode Reader (RS232C)	1 Serial Splitter Cable (RS232C)	23.6 cm x 10.8 cm x 7.8 cm; 0.482 kg 9.3" x 4.3" x 3.1"; 17.0 oz	22
Printer Paper Rolls	U121-101	4 Printer Paper Rolls	Thermal Paper (0.06 m x 20 m): 200 results/roll Sticker Paper (0.06 m x 9 m): 100 results/roll	12.0 cm x 12.0 cm x 6.5 cm; 0.360 kg 4.7" x 4.7" x 2.6"; 12.7oz 63.0 cm x 37.0 cm x 30.0 cm; 19.4 kg 24.8" x 14.6" x 11.8"; 684.3 oz 12.0 cm x 12.0 cm x 6.5 cm; 0.40 kg 4.7" x 4.7" x 2.6"; 14.1oz 63.0 cm x 37.0 cm x 30.0 cm; 21.4 kg 24.8" x 14.6" x 11.8"; 684.3 oz; 754.9 oz	50
U500 Data Transfer Kit	U221-131 <sup>✓</sup>	1 Data Transfer Cable (RS232C)	1 Package Insert	16.0 cm x 13.0 cm x 3.5 cm; 0.147kg 6.3" x 5.1" x 1.4"; 5.2 oz	8

**We also offer other rapid diagnostic and medical products:**

Blood Glucose Monitoring Systems, Immunoassay EIA/ELISA and more.

✓ CE Marked for sale in the European Community **CE**

† Cleared for US 510(k)

## Nutrient Agar

M001

### Intended use

Nutrient Agar is used as a general purpose medium for the cultivation of less fastidious microorganisms, can be enriched with blood or other biological fluids.

### Composition\*\*

Ingredients	Gms / Litre
Peptone	5.000
Sodium chloride	5.000
HM peptone B <sup>#</sup>	1.500
Yeast extract	1.500
Agar	15.000
Final pH ( at 25°C)	7.4±0.2

\*\*Formula adjusted, standardized to suit performance parameters

# - Equivalent to Beef extract

### Directions

Suspend 28.0 grams in 1000 ml purified / distilled water. Heat to boiling to dissolve the medium completely. Sterilize by autoclaving at 15 lbs pressure (121°C) for 15 minutes. Cool to 45-50°C. If desired ,the medium can be enriched with 5-10% blood or other biological fluids. Mix well and pour into sterile Petri plates.

### Principle And Interpretation

Nutrient media are basic culture media used for maintaining microorganisms, cultivating fastidious organisms by enriching with serum or blood and are also used for purity checking prior to biochemical or serological testing (1,2). Nutrient Agar is ideal for demonstration and teaching purposes where a more prolonged survival of cultures at ambient temperature is often required without risk of overgrowth that can occur with more nutritious substrate. This relatively simple formula has been retained and is still widely used in the microbiological examination of variety of materials and is also recommended by standard methods. It is one of the several non-selective media useful in routine cultivation of microorganisms (3,4). It can be used for the cultivation and enumeration of bacteria which are not particularly fastidious. Addition of different biological fluids such as horse or sheep blood, serum, egg yolk etc. makes it suitable for the cultivation of related fastidious organisms. Peptone, HM peptone B and yeast extract provide the necessary nitrogen compounds, carbon, vitamins and also some trace ingredients necessary for the growth of bacteria. Sodium chloride maintains the osmotic equilibrium of the medium.

### Type of specimen

Clinical samples - faeces, urine ; Food and dairy samples; Water samples

### Specimen Collection and Handling:

For clinical samples follow appropriate techniques for handling specimens as per established guidelines (5,6).

For food and dairy samples, follow appropriate techniques for sample collection and processing as per guidelines (3,4,7).

For water samples, follow appropriate techniques for sample collection, processing as per guidelines and local standards (8).

After use, contaminated materials must be sterilized by autoclaving before discarding.

### Warning and Precautions :

In Vitro diagnostic use. For professional use only. Read the label before opening the container. Wear protective gloves/protective clothing/eye protection/ face protection. Follow good microbiological lab practices while handling specimens and culture. Standard precautions as per established guidelines should be followed while handling clinical specimens. Safety guidelines may be referred in individual safety data sheets.

### Limitations :

1. Individual organisms differ in their growth requirement and may show variable growth patterns on the medium.

2. Each lot of the medium has been tested for the organisms specified on the COA. It is recommended to users to validate the medium for any specific microorganism other than mentioned in the COA based on the user's unique requirement.

## Performance and Evaluation

Performance of the medium is expected when used as per the direction on the label within the expiry period when stored at recommended temperature.

### Quality Control

#### Appearance

Cream to yellow homogeneous free flowing powder

#### Gelling

Firm, comparable with 1.5% Agar gel

#### Colour and Clarity of prepared medium

Light yellow coloured clear to slightly opalescent gel forms in Petri plates

#### Reaction

Reaction of 2.8% w/v aqueous solution at 25°C. pH : 7.4±0.2

#### pH

7.20-7.60

#### Cultural Response

Cultural characteristics observed after an incubation at 35-37°C for 18-48 hours.

Organism	Inoculum (CFU)	Growth	Recovery
<i>Escherichia coli</i> ATCC 25922 (00013*)	50-100	good-luxuriant	≥70%
<i>Pseudomonas aeruginosa</i> ATCC 27853 (00025*)	50-100	good-luxuriant	≥70%
<i>Salmonella</i> Typhi ATCC 6539	50-100	good-luxuriant	≥70%
<i>Staphylococcus aureus</i> subsp. <i>aureus</i> ATCC 25923 (00034*)	50-100	good-luxuriant	≥70%
<i>Streptococcus pyogenes</i> ATCC 19615	50-100	good-luxuriant	≥70%
<i>Salmonella</i> Enteritidis ATCC 13076 (00030*)	50-100	good-luxuriant	≥70%
<i>Salmonella</i> Typhimurium ATCC 14028 (00031*)	50-100	good-luxuriant	≥70%
<i>Yersinia enterocolitica</i> ATCC 9610 (00038*)	50-100	good-luxuriant	≥70%
<i>Yersinia enterocolitica</i> ATCC 23715 (00160*)	50-100	good-luxuriant	≥70%

Key : (\*) Corresponding WDCM numbers.

### Storage and Shelf Life

Store between 10-30°C in a tightly closed container and the prepared medium at 20-30°C. Use before expiry date on the label. On opening, product should be properly stored dry, after tightly capping the bottle in order to prevent lump formation due to the hygroscopic nature of the product. Improper storage of the product may lead to lump formation. Store in dry ventilated area protected from extremes of temperature and sources of ignition. Seal the container tightly after use. Product performance is best if used within stated expiry period.

### Disposal

User must ensure safe disposal by autoclaving and/or incineration of used or unusable preparations of this product. Follow established laboratory procedures in disposing of infectious materials and material that comes into contact with clinical sample must be decontaminated and disposed of in accordance with current laboratory techniques (5,6).

## Reference

- 1.Lapage S., Shelton J. and Mitchell T., 1970, Methods in Microbiology', Norris J. and Ribbons D., (Eds.), Vol. 3A, Academic Press, London.
- 2.MacFaddin J. F., 2000, Biochemical Tests for Identification of Medical Bacteria, 3rd Ed., Lippincott, Williams and Wilkins,Baltimore.
- 3.American Public Health Association, Standard Methods for the Examination of Dairy Products, 1978, 14th Ed., Washington D.C.
- 4.Salfinger Y., and Tortorello M.L., 2015, Compendium of Methods for the Microbiological Examination of Foods, 5th Ed., American Public Health Association, Washington, D.C.
- 5.Isenberg, H.D. Clinical Microbiology Procedures Handbook 2nd Edition.
- 6.Jorgensen, J.H., Pfaller, M.A., Carroll, K.C., Funke, G., Landry, M.L., Richter, S.S and Warnock., D.W. (2015) Manual of Clinical Microbiology, 11th Edition. Vol. 1.
- 7.Wehr H. M. and Frank J. H., 2004, Standard Methods for the Microbiological Examination of Dairy Products, 17th Ed.,APHA Inc., Washington, D.C.
- 8.Baird R.B., Eaton A.D., and Rice E.W., (Eds.), 2015, Standard Methods for the Examination of Water and Wastewater, 23rd ed., APHA, Washington, D.C.


Revision : 06/2022

HiMedia Laboratories Pvt. Limited,  
Plot No.C-40, Road No.21Y,  
MIDC, Wagle Industrial Area,  
Thane (W) -400604, MS, India

**IVD** *In vitro* diagnostic  
medical device

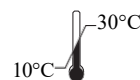**Storage temperature**

CEpartner4U, Esdoornlaan 13,  
3951DB Maarn, NL  
www.cepartner4u.eu

**CE Marking**

**Do not use if  
package is damaged**

### Disclaimer :

User must ensure suitability of the product(s) in their application prior to use. Products conform solely to the information contained in this and other related HiMedia™ publications. The information contained in this publication is based on our research and development work and is to the best of our knowledge true and accurate. HiMedia™ Laboratories Pvt Ltd reserves the right to make changes to specifications and information related to the products at any time. Products are not intended for human or animal or therapeutic use but for laboratory,diagnostic, research or further manufacturing use only, unless otherwise specified. Statements contained herein should not be considered as a warranty of any kind, expressed or implied, and no liability is accepted for infringement of any patents.

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
accuracy and quality as a science

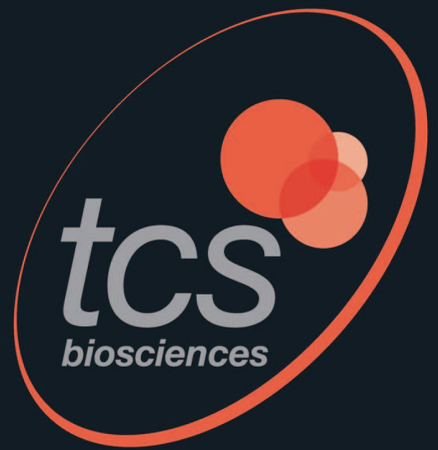Selectrol®  
Technical  
Guide

**Selectrol®** : Manufactured under licence from Public Health England Culture Collections

## SELECTROL® - FREEZE-DRIED ORGANISMS IN A DISC

**Quality control of microbial characterisation tests, culture media and antimicrobial susceptibility determinations is best accomplished by the use of microorganisms with well-documented and stable phenotypic and genotypic characteristics.**

Bacterial and fungal strains have been selected and recommended by expert bodies, such as **EUCAST**, **CLSI** and the European Pharmacopoeia, on the basis of their suitability for monitoring test performance and ensuring the validity of results for testing used in clinical, food, pharmaceutical, water and veterinary laboratories.

Products derived from the cultures in the collections should be manufactured using the minimum number of sub-cultures, to minimise the possibility of alterations to the phenotype due to mutations. See also page 14.

Selectrol strains are manufactured exclusively from Public Health England Culture Collections (NCTC® and NCPF®) and are first generation subcultures, unlike many products on the market which are 2<sup>nd</sup>, 3<sup>rd</sup> or 4<sup>th</sup> generation subcultures. They are preserved by long-term storage as freeze-dried cells in order to minimise any alterations to the phenotype caused by mutations.

### Passages

A Selectrol® disc is a first generation subculture from a **master culture** sourced from Public Health England Culture Collections, and is designed to be used to obtain **working stock** cultures for use in testing. It is generally accepted that no more than a total of five passages should be made from the **master culture**, in order to avoid genetic drift and mutant selection. Therefore, no more than four passages (fresh cultures) from the **working stock** should be made.

### Shelf life

For most strains, Selectrol® discs are guaranteed to contain at least 10<sup>6</sup> organisms at the time of purchase; this number is sufficient to ensure that when the discs are used and stored as directed there will be viable organisms cultivable up to the stated end of the shelf life, which is usually 9 months from the time the vial is first opened.

### Quality Control

Selectrol® batches are tested in our UKAS accredited testing laboratory number 2496. A test report for each batch of Selectrol® can be accessed via our website. The reporting of Selectrol® test results via the website comes under our UKAS accreditation.

Selectrol® cultures are rigorously tested to confirm identity, to confirm the possession of essential phenotypic characteristics and to exclude contamination with other organisms. Photographic evidence of the test results is retained for each batch, along with retained appropriately stored samples.

## Glossary

**AMRHAI:** Antimicrobial Resistance and Healthcare Associated Infections reference unit

**ATCC®:** American Type Culture Collection. ATCC® strains are listed for reference only. ATCC® is a registered trademark of the American Type Culture Collection.

**BSAC:** British Society for Antimicrobial Chemotherapy - Now superseded by EUCAST

**CLSI:** Clinical Laboratory Standards Institute. (USA)

**CPE:** Carbapenemase Producing Enterobacteriaceae

**CRE:** Carbapenem Resistant Enterobacteriaceae

**Culture collection:** Cultures of fully characterised organisms maintained in such a way as to minimise sub-culturing. See page 14.

**ESBL:** Extended Spectrum Beta-Lactamase-producing organism.

**EUCAST:** European Committee on Antimicrobial Susceptibility Testing.

**First generation derivative:** A single passage from a master culture, for example a Selectrol® disc.

**Master culture:** Culture derived from a reference culture vial.

**NCPF®:** National Collection of Pathogenic Fungi. NCPF® is a registered trademark of Public Health England.

**NCTC®:** National Collection of Type Cultures. NCTC® is a registered trademark of Public Health England.

**Passage:** An equivalent term for a subculture.

**PHE:** Public Health England.

**Reference cultures:** Quality control strains selected on the basis of their phenotypic biochemical and antimicrobial susceptibility characteristics to be used as controls in microbiological testing. These are obtained as freeze-dried vials from culture collections.

**Stock culture:** Cultures derived from a Selectrol® disc, which can be stored for up to a week, usually on agar slants.

**Working cultures:** Stock cultures further sub-cultured to provide 18-24 hour growth for use in testing.

**WDCM:** World Data Centre for Microorganisms

**WFCC:** World Federation for Culture Collections

## SIGNIFICANT PROPERTIES AND USES OF SELECTROL® ORGANISMS

### ***Aspergillus brasiliensis*** (formerly *Aspergillus niger*):

MM94 – NCPF® 2275 / ATCC® 16404 / WDCM 00053 – used in pharmaceutical industry for testing media and preservatives. Colonies are initially white or yellowish and on the reverse greyish or greenish-yellow. Sporing heads on the colony surface are initially pale, becoming dark brown to black. Sporulation may be inhibited in sealed plates.

### ***Bacillus cereus***:

MM21 – NCTC® 10320 / ATCC® 9634 / WDCM 00001 (recently renamed *Bacillus toyonensis*) – ISO 11133 recommended media and ID test control organism.

MM86 – NCTC® 7464 / ATCC® 10876 – **PHE** recommended media and ID test control organism.

### ***Bacillus subtilis*** (*Bacillus subtilis* subsp. *spizizenii*):

MM29 – NCTC® 10400 / ATCC® 6633 / WDCM 00003 – used in antibiotic assays (fully sensitive), **PHE** recommended media and ID test control organism.

### ***Bacteroides fragilis***:

MM44 – NCTC® 9343 / ATCC® 25285 – type strain, **PHE** recommended strain for media and sensitivity test control.

### ***Campylobacter jejuni*** (*Campylobacter jejuni* subsp. *jejuni*):

MM82 – NCTC® 11322 / ATCC® 29428 / WDCM 00156 – **PHE** recommended strain for media control.

MM36 – NCTC® 11351 / ATCC® 33560 – **EUCAST** recommended strain for susceptibility testing.

### ***Candida albicans***:

MM28 – NCPF® 3255 / ATCC® 2091 / WDCM 00055 – sensitivity control / industrial use.

MM42 – NCPF® 3179 / ATCC® 10231 / WDCM 00054 – pharmaceutical / media testing / **PHE** recommended strain for media control.

### **CRE ≡ ‘Carbapenem Resistant Enterobacteriaceae’ / CPE ≡ ‘Carbapenemase Producing Enterobacteriaceae’**

There are 5 carbapenemases which are currently a significant problem in the UK – KPC, OXA-48, IMP, NDM and VIM – and PHE recommend that all clinically-significant Gram-negative bacteria should be routinely screened for carbapenemase production, using a recommended carbapenem<sup>2</sup> such as ertapenem or meropenem. Resistant isolates may be investigated further to determine which resistance mechanism is involved using the Modified Hodge Test, MALDI-TOF, PCR or a reference laboratory.

MM55 *Klebsiella pneumoniae* – NCTC® 13440 – produces a Class B VIM-1 Carbapenemase.

MM56 *Klebsiella pneumoniae* – NCTC® 13443 – produces a Class B NDM-1 Carbapenemase.

MM58 *Klebsiella pneumoniae* – NCTC® 13438 – produces a Class A KPC-3 Carbapenemase.

MM59 *Klebsiella pneumoniae* – NCTC® 13442 – produces a Class D OXA-48 Carbapenemase.

MM57 *Escherichia coli* – NCTC® 13476 – produces a Class B IMP Carbapenemase.

MM33 *Escherichia coli* – NCTC® 10418 / ATCC® 10536 – recommended by **PHE** as a negative control for CRE testing.

***Citrobacter freundii:***

MM27 – NCTC® 9750 / ATCC® 8090 – type strain.

***Clostridium perfringens:***

MM45 – NCTC® 8237 / ATCC® 13124 / WDCM 00007 – type strain. **PHE** recommended strain for food testing (Tryptose Sulphite Cycloserine agar – lactose and gelatin positive) and sensitivity test control. *Clostridium perfringens* is listed in Schedule 5 of the Anti-terrorism, Crime and Security Act 2001, and should be securely stored in accordance with the guidelines of the Act. However, MM45 is a type A strain, which does not produce the lethal epsilon toxin of potential interest to bioterrorists.

***Clostridium sporogenes:***

MM31 – NCTC® 532 / ATCC® 19404 / WDCM 00008 – used for media control. **PHE** recommended strain for media QC (lactose gelatin medium for ID of *C. perfringens* lactose negative and gelatin positive).

***Enterobacter aerogenes:***

MM26 – NCTC® 10006 / ATCC® 13048 / WDCM 00175 – type strain; used in water, paint and adhesive testing.

***Enterobacter cloacae:***

MM01 – NCTC® 13380 / ATCC® 23355 / WDCM 00082 – disinfectant control, media testing.

MM51 – NCTC® 13406 – **PHE** recommended strain for QC of AmpC (de-repressed) detection.

***Enterococcus faecalis:***

MM52 – NCTC® 13379 / ATCC® 51299 / WDCM 00085 – is vancomycin resistant (low-level VanB mediated) and also shows high-level resistance to aminoglycosides. It is used to confirm methodologies used to detect these resistances are working correctly. Lancefield group D.

MM17 – NCTC® 775 / ATCC® 19433 / WDCM 00009 – used in water industry and QC. **PHE** recommended strain for media control. Fully sensitive. Lancefield group D.

MM18 – NCTC® 12697 / ATCC® 29212 / WDCM 00087 – is fully sensitive to vancomycin and gentamicin. **PHE** recommended positive control strain for aesculin test. **CLSI, EUCAST** recommended media control for sulpha / trimethoprim testing and general susceptibility testing control. Lancefield group D.


***Enterococcus hirae:***

MM35 – NCTC® 13383 / ATCC® 10541 / WDCM 00011 – disinfectant control. Used in microbiological assays. Colonies are alpha-haemolytic on sheep blood agar.

***Escherichia coli* strains:**

MM02 – NCTC® 12241 / ATCC® 25922 / WDCM 00013 – **EUCAST**, **CLSI**, **PHE** recommended control strain for susceptibility testing (fully sensitive). Exhibits 2 colony types – the most prevalent type is slightly irregular, smooth and translucent. The secondary type appears more opaque. It is preferable to maintain cultures on agar as passage in broth can result in a change in MIC levels.

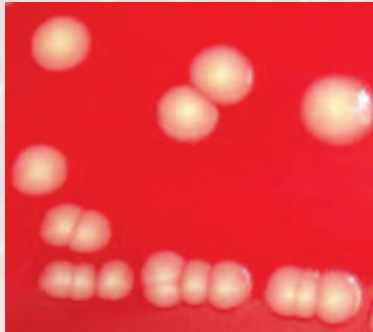MM57 – NCTC® 13476 – CRE testing control; produces a Class B IMP Carbapenemase.

MM33 – NCTC® 10418 / ATCC® 10536 – (**PHE** recommended alternative to NCTC 12241) fully sensitive control strain. **PHE** recommended positive control for indole test, ONPG test, negative control for oxidase test, **PHE** recommended negative control for CRE and ESBL testing.

MM24 – NCTC® 11954 / ATCC® 35218 – beta-lactamase positive strain. **CLSI** recommended strain for susceptibility testing ONLY for penicillin / beta-lactamase inhibitor combinations. Sensitive to amoxicillin / clavulanic acid.

MM75 – NCTC® 9001 / ATCC® 11775 / WDCM 00090 – used in water / chemical industry. **PHE** recommended strain for media QC.

MM93 – NCTC® 12900 / ATCC® 700728 / WDCM 00014 – O157 strain (non-toxigenic). **PHE** recommended strain for media QC.

MM63 – NCTC® 11560 – beta-lactamase positive strain.

MM38 – NCTC® 12923 / ATCC® 8739 / WDCM 00012 – used in pharmaceutical / water industry. Three colony types: A) Entire, glistening, smooth and translucent. B) Entire, glistening smooth and opaque. C) Irregular, rough and translucent. The rough colonies appear after 48 hours incubation.

MM34 – NCTC® 13846 – Possesses the plasmid-mediated *mcr-1* colistin resistance mechanism gene and is recommended by **PHE** and **EUCAST** as a control for tests to detect this increasingly prevalent resistance, in conjunction with NCTC® 12241 / ATCC® 25922 (Selectrol strain MM02) as a negative control.


***Haemophilus influenzae* strains:**

MM81 - NCTC® 12699 / ATCC® 49247 – is a ‘BLNAR’ strain – (beta-lactamase non-producing ampicillin / amoxycillin resistant). These strains are important clinically because the susceptibility results obtained using conventional testing procedures maybe misleading in the case cephalosporins. **PHE**, **CLSI** recommended QC strain for susceptibility testing media.

MM98 – NCTC® 11931 – a fully sensitive strain. **PHE** recommended strain for porphyrin synthesis test, chocolate agar control.

MM100 – NCTC® 8468 / ATCC® 9334 / CCUG 23946 – another fully sensitive strain, which reportedly gives results which are easier to interpret when Mueller-Hinton medium is used in preference to Iso-Sensitest medium. MIC for amoxycillin is 0.5 mg/l.

MM37 – NCTC® 12975 / ATCC® 49766 – recommended by **EUCAST**.

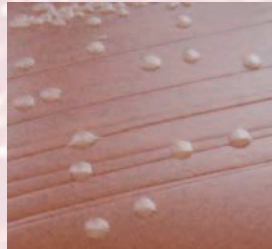
***Klebsiella* strains:**

MM04 *Klebsiella pneumoniae* – NCTC® 9633 / ATCC® 13883 / WDCM 00097 – type strain. Two colony types may be seen. The predominant type is entire and opaque. The secondary type is slightly smaller and translucent.

MM83 *Klebsiella pneumoniae* – NCTC® 13368 / ATCC® 700603 – ESBL-producing strain used as control for ESBL testing. There are two colony types.

MM55 *Klebsiella pneumoniae* – NCTC® 13440 – CRE testing control; produces a Class B VIM-1 Carbapenemase.

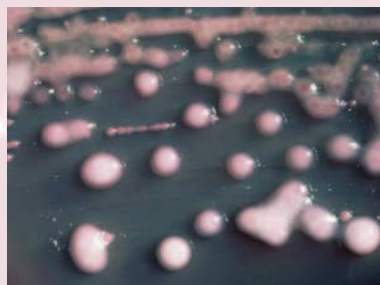MM56 *Klebsiella pneumoniae* – NCTC® 13443 – CRE testing control; produces a Class B NDM-1 Carbapenemase.

MM58 *Klebsiella pneumoniae* – NCTC® 13438 – CRE testing control; produces a Class A KPC-3 Carbapenemase.

MM59 *Klebsiella pneumoniae* – NCTC® 13442 – CRE testing control; produces a Class D OXA-48 Carbapenemase.

MM88 *Klebsiella aerogenes* (*Raoultella planticola*) – NCTC® 9528 – used in water / pharmaceutical industry. **PHE** recommended negative control for Tryptone Bile X-Glucuronide agar and Yeast Extract agar.

***Lactobacillus brevis:***

MM76 – NCTC® 13386 / ATCC® 8287 – used in food industry.

***Legionella pneumophila serogroup 1:***

MM08 – NCTC® 11192 / ATCC® 33152 / WDCM 00107 – derived from strain isolated from first recognised outbreak of legionellosis in Philadelphia at the Legionnaires' Convention 1976

***Listeria innocua:***

MM92 – NCTC® 11288 / ATCC® 33090 / WDCM 00017 – type strain. Non-pathogenic.

***Listeria monocytogenes:***

MM87 – NCTC® 11994 / WDCM 00019 – type strain, **PHE** recommended positive control strain for Listeria detection in food. Serotype 4b, most common serovar isolated from human infections.

MM48 – NCTC® 7973 / ATCC® 35152 / WDCM 00109 – produces 2 phenotypes, one is beta-haemolytic and virulent, the other non-haemolytic and non-virulent. Serovar 1/2a.

MM77 – NCTC® 13372 / ATCC® 7644 – used in food microbiology Q.C. Colonies exhibit beta-haemolysis on sheep blood agar.

***Neisseria gonorrhoeae:***

MM96 – NCTC® 12700 / ATCC® 49226 – has low-level, but clinically relevant, resistance to penicillin – MIC of penicillin is 0.5 mg/l. **PHE** recommended control for susceptibility testing – methodology assesses the ability of testing to detect resistance rather than sensitivity; this strain has low-level, but clinically relevant, resistance to penicillin – MIC of penicillin is 0.5 mg/l. Some variation in size and texture of colonies may be observed. Increased CO<sub>2</sub> is helpful in growth.

MM05 – NCTC® 8375 / ATCC® 19424 – is fully sensitive – MIC of penicillin is 0.06 mg/l. **PHE** recommended strain for media QC.

***Proteus mirabilis:***

MM43 – NCTC® 13376 / ATCC® 14153 – pharmaceutical / disinfectant / media control.

MM68 – NCTC® 10975 – media control. **PHE** recommended control for motility test.


#### ***Proteus vulgaris:***

MM09 – NCTC® 4175 / ATCC® 13315 – was the type strain, but is atypical and has been recognised as a separate species – *Proteus hauseri* – it is used for media control. Colonies are glistening with spreading edges.

#### ***Pseudomonas aeruginosa* strains:**

MM10 – NCTC® 12903 / ATCC® 27853 / WDCM 00025 – is fully sensitive to anti-pseudomonal antibiotics (EUCAST susceptibility test control). 2 colony types may be observed: A) predominantly flat, spreading edges and rough surface; B) small and compact. Produces both fluorescein and pyocyanin pigments.

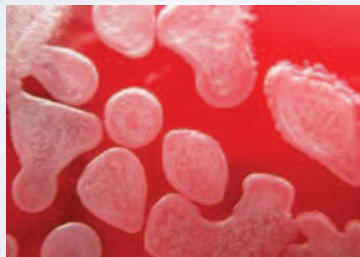MM65 – NCTC® 10662 / ATCC® 25668 / WDCM 00114 – is fully sensitive. PHE recommended control strain for media control

MM40 – NCTC® 12924 / ATCC® 9027 / WDCM 00026 – used in water industry / disinfectant testing. Colonies on agar plates are entire, glistening and mucoid with a grainy surface. This strain also produces both fluorescein and pyocyanin pigments.

MM41 – NCTC® 13359 / ATCC® 15442 – used in water industry / disinfectant testing. May produce up to 3 different colony types. Pyocyanin is not produced.

#### ***Rhodococcus equi:***

MM97 – NCTC® 1621 / ATCC® 6939 / WDCM 00028 – type strain.

#### ***Saccharomyces cerevisiae:***

MM73 – NCPF® 3178 – PHE recommended strain for food testing and enumeration of yeasts and moulds.

MM50 — NCTC® 10716 / WDCM 00058 – used for QC of culture media and for antifungal susceptibility testing.

#### ***Salmonella* serotypes:**

MM11 *Salmonella* Typhimurium – NCTC® 12023 / ATCC® 14028 / WDCM 00031 – (1,4,5,12: i: 1,2) Used for media/test QC. This is a common serotype from animals and from human infections.

The strains listed below are unusual serotypes, used to avoid any chance of confusion with strains commonly found in animals, food, etc, and are used to control media and detection methods in the food industry:

MM89 *Salmonella* Poona – NCTC® 4840 – (13,22: z: 1,6) PHE recommended control strain for food testing.

MM84 *Salmonella* Nottingham – NCTC® 7832 – (16: d: e,n,z15) PHE recommended control for water testing.

#### ***Serratia marcescens:***

MM12 – NCTC® 13382 / ATCC® 8100 – used for disinfectant testing. PHE recommended negative control for indole test. Colonies are entire, glistening, smooth and translucent. Non-pigmented.

### ***Staphylococcus aureus:***

#### (A) Fully sensitive:

MM85 – NCTC® 6571 / ATCC® 9144 / WDCM 00035 – historically used for susceptibility testing ('Oxford staph'), but largely superseded by MM13 as it has unusually low MIC's and so is unrepresentative of normal range of Staph aureus strains. Sensitive to penicillin and ceftazidime / methicillin / oxacillin. **PHE** recommended coagulase, DNase and catalase positive control.

MM13 – NCTC® 12981 / ATCC® 25923 / WDCM 00034 – used in susceptibility and media testing/QC. Fully sensitive to all anti-staphylococcal antibiotics (including penicillin and methicillin / oxacillin). It is preferable to maintain cultures on agar as passage in broth can result in a change in MIC levels. Colonies are circular white to cream, convex to flat in elevation. After 48 hours incubation a few grey/translucent variants may be noted. Beta-haemolytic on sheep blood agar.

#### B) Penicillin resistant:

MM14 – NCTC® 12973 / ATCC® 29213 / WDCM 00131 – used for susceptibility testing, especially for automated methodology. **EUCAST**, **CLSI** strain. Sensitive to ceftazidime / methicillin / oxacillin. Penicillin resistant – weak beta-lactamase producer. Colonies are beta-haemolytic, and a golden-orange colour.

MM30 – NCTC® 7447 / ATCC® 6538P / WDCM 00033 – used for susceptibility testing/antibiotic assay, disinfectant testing. Ceftazidime / methicillin / oxacillin sensitive. Penicillin resistant. Colonies are weakly beta-haemolytic, coagulase positive and beta-lactamase negative.

#### (C) MRSA (ceftazidime / methicillin / oxacillin resistant):


MM91 – NCTC® 13373 / ATCC® 43300 / WDCM 00211 (MRSA) – Possesses mecA gene but is hetero-resistant, (so as few as one per thousand cells demonstrate the resistance) and consequently has low-level ceftazidime / oxacillin / methicillin resistance (4.0 mg/l MIC of oxacillin, 8.0 mg/l MIC of ceftazidime – methicillin sensitive strains have MIC of 0.12-0.5 for oxacillin and 1-4 for ceftazidime.); it is used to confirm testing procedures for methicillin resistance are working and provides a more stringent test than testing with an MRSA which shows homogeneous resistance and has a much higher MIC. This organism will have a zone of inhibition reduced in size compared to a fully ceftazidime / oxacillin / methicillin sensitive strain (such as MM13). **CLSI** recommended strain for MRSA testing. There are two colony types: 1) Beta-haemolytic with a slight yellow tint. 2) Non-haemolytic and white.

MM64 – NCTC® 12493 / WDCM 00212 (MRSA) – possesses mecA gene and shows homogeneous resistance with MIC of >64 for methicillin, which produces high-level ceftazidime / methicillin / oxacillin resistance. **EUCAST** recommended strain. Instances have been reported where loss of the mecA gene has occurred during storage.

#### D) Other:

MM46 – NCTC® 10788 / ATCC® 6538 / WDCM 00032 – used in pharmaceutical industry for testing disinfectants etc. Usually yellow pigmented colonies, or can produce a white colonial variant. Beta-haemolytic.

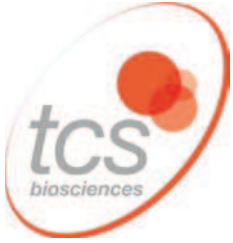***Staphylococcus epidermidis:***

MM15 – NCTC® 13360 / ATCC® 12228 / WDCM 00036 – used for media control / antibiotic assay. Colonies are small and beta-haemolytic.

***Streptococcus agalactiae:*** (Beta-haemolytic Streptococcus group B)

MM16 – NCTC® 8181 / ATCC® 13813 – type strain, used for QC. PHE recommended negative control for aesculin test.

***Streptococcus pneumoniae*** strains:

MM95 – NCTC® 12977 / ATCC® 49619 – has low-level, but clinically relevant, resistance to penicillin – this organism is used to assess detection of resistance rather than sensitivity. PHE recommended positive control for bile solubility test. CLSI, EUCAST recommended control strain for susceptibility testing. Serotype 19F.

MM19 – NCTC® 12695 / ATCC® 6303 – is fully sensitive. Colonies are mucoid and alpha-haemolytic. A few colonies may have an irregular edge. Serotype 3.

***Streptococcus pyogenes:***

MM20 – NCTC® 12696 / ATCC® 19615 – used for QC and media testing. Lancefield group A, beta-haemolytic. PHE recommended blood agar control.

***Vibrio parahaemolyticus:***

MM06 – NCTC® 10885 / WDCM 00185 – used for QC of media and ID testing. PHE recommended strain used mainly in the food industry.

***Yersinia enterocolitica:***

MM80 – NCTC® 12982 / ATCC® 9610 / WDCM 00038 – type strain, used for media control. Serotype O:8, which is a pathogenic serotype, commonest in USA.

**References:**


- 1 European Committee on Antimicrobial Susceptibility Testing (EUCAST). Routine and Extended Internal Quality Control for MIC Determination and Disc Diffusion. Version 7.0 - 01.01.2017.
- 2 UK Standards for Microbiology Investigations. Example Reference Strains for Microbiology Investigations Test Procedures: Bacteriology—Test Procedures | TP 1 | Issue No. 2 | 05.01.2015. Public Health England (PHE).
- 3 Performance Standards for Antimicrobial Disc Susceptibility Tests: Approved Standard—11th Edition. Clinical and Laboratory Standards Institute (CLSI).

## How to use Selectrol®


Always warm the vial to ambient temperature before opening.

Be sure to use non-selective culture media to revive the organisms.

For the more fastidious organisms, such as anaerobes, it is generally better to use agar rather than broth for revival.

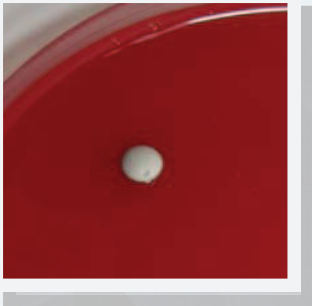
Place disc on suitable growth medium such as blood agar

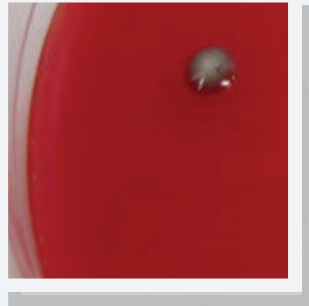
Leave disc for a few minutes to liquefy, then spread plate and incubate to produce isolated colonies

Obtain a stock culture which can be used to prepare an inoculum for biochemical and antibiotic susceptibility tests

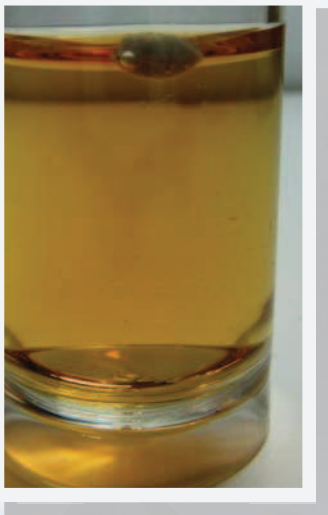Place disc in a small volume of a suitable broth medium such as brain-heart infusion

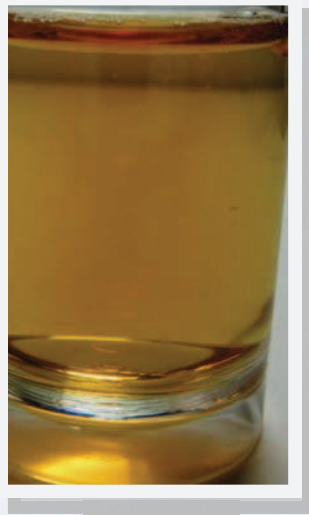Allow disc a few minutes to dissolve, then spread aliquot onto a plate of suitable growth medium

## Out-of-specification results


Laboratories use Selectrol® for Quality Control of culture media, biochemical identification tests and antimicrobial susceptibility testing. When a laboratory test result, an MIC or biochemical reaction, is unexpected or out-of-specification, the test should first be repeated to confirm it; an out-of-specification result is an indication that the testing procedure should be reviewed; it is not, in the first instance, a sign of a problem with the control organism.

If incorrect results are obtained on retesting, the explanation could be:

- The test procedure was not followed correctly – check standard operating procedures
- There is an instrumentation error – check calibration, mechanical functioning, etc
- There is a problem with the consumables – out of date, incorrect storage, etc
- The culture of the control organism has become contaminated

### Technical Support

If no explanation for out-of-spec results can be found, but repeated tests still give unacceptable results, please contact TCS and / or your relevant reference laboratory or instrument manufacturer for advice. For example, contact AMRHAI at Colindale, London if MIC results are consistently outside the acceptable range. Please retain any remaining discs of organisms about which you have concerns so they can be returned to TCS and investigated alongside retained samples.
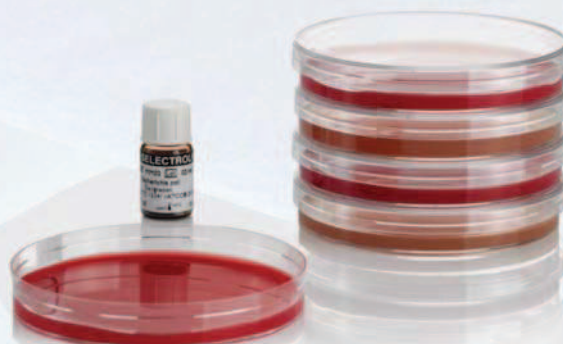
## Preparing QC and Validation Spikes from Selectrol®

### Preparing the spike

- Place a Selectrol® disc in Brain Heart Infusion (BHI) broth\* or equivalent, and culture (typically for 18 hours) at the appropriate temperature for the organism (typically 37°C)
- Assume the count in the broth to be  $10^8$  organisms per ml ----- (A)
- Mix and transfer 100  $\mu$ l of (A) to 100 ml of saline or  $\frac{1}{4}$  strength Ringer's solution -- (B)
- Mix and transfer 100  $\mu$ l of (B) to 10 ml of saline or  $\frac{1}{4}$  strength Ringer's solution --- (C)
- Mix and transfer 100  $\mu$ l of (C) to your homogenised food sample.

### Verifying the inoculum

- Pipette 5 x 10  $\mu$ l drops from (C) onto each of two agar plates for Miles and Misra counts.

### Using the assumptions and dilutions above:

- (A) contains  $10^8$  organisms per ml
- (B) contains  $10^5$  organisms per ml
- (C) contains  $10^3$  organisms per ml

### If the Miles and Misra counts indicate that the required count was not achieved:

- If the count was too high by a factor of 10, reduce the volume transferred from (A) to (B) from 100  $\mu$ l to 10  $\mu$ l
- If the count was too low by a factor of 10, increase the volume transferred from (A) to (B) from 100  $\mu$ l to 1 ml.

Keep a record of the correct dilutions for each organism type for future use. You will find that this method is very repeatable.

\*Note: BHI broth will work for most of the Selectrol® organisms; however, for fastidious organisms an appropriate culture broth must be selected, e.g. Fastidious Anaerobe Broth for strictly anaerobic organisms.

## Culture Collections

Cultures of microorganisms have been deposited and subsequently maintained in 589 collections in 68 countries, and many of the cultures are derived from the same original isolate; the history of each organism, its properties and names of the culture collections which hold it are detailed in the relevant catalogues and websites.

Some of the organisms have been selected and recommended by expert organisations to be supplied as controls for microbiological tests, and when the identical cultures are present in more than one collection they will have a specific designation for each, incorporating the abbreviation for the collection and a reference number.

For example:- *Staphylococcus aureus* NCTC 7447, widely recommended as a control for antimicrobial susceptibility testing, is held in 30 collections, and consequently the phenotypically and genotypically identical organism has 30 different references, such as ATCC 6538P, CIP 53.156, DSM 346 and so on.

In an effort to minimise potential confusion and help users find local sources of reference strains, the WFCC and the WDCM initiated a system that ascribes each recommended QC strain a reference number (WDCM 00001 onwards), cites all collections that contain it and provides contact details and each collection's unique reference. For example, the strain of *Staphylococcus aureus* NCTC 7447 (Selectrol® strain MM33) mentioned above is designated WDCM 00033.

### ***Staphylococcus aureus* WDCM 00033**

AHU 1142; **ATCC™ 6538P**; BCRC 10451 ; BTCC 209P; BU 395; CCM 2022; CCTM 596; CCUG 1828; CECT 240; CIP 53.156; CN 3784; CNCTC Mau 28/58; DSM 346; FIRDI 451; IAM 1011; IAM 12082; IEM Mau 28/58; IFO 12732; IFO 3061; IID 671; IMET 10904; JCM 2151; LMG 8195; NCIMB 8625; **NCTC 7447**; NRRL B-313; OUT 8232; PCI 1209; PZH 8/54; RIMD 3109007; VNIIA 209P;

Products derived from the cultures in the collections should be manufactured using the minimum number of sub-cultures, to minimise the possibility of alterations to the phenotype due to mutations. Ideally, as in the case of **Selectrol®**, a single sub-culture only is used, so the **Selectrol®** product is a 'first generation derivative' of a culture supplied by NCTC, and will be identical with regard to its properties and suitability for use in QC applications to a culture of the particular organism obtained from any of the other WDCM listed culture collections.

---

Every effort has been made to ensure the accuracy of the information in this document, however TCS makes no warranties, expressed or implied, regarding errors or omissions and assumes no legal liability or responsibility for loss or damage resulting from the use of information contained within.

## Selectrol Strain Index

Strain Name	Designation	Code	WDCM
<i>Aspergillus brasiliensis</i>	NCPF <sup>®</sup> 2275 / ATCC <sup>®</sup> 16404	MM94	00053
<i>Bacillus cereus</i>	NCTC <sup>®</sup> 10320 / ATCC <sup>®</sup> 9634	MM21	00001
<i>Bacillus cereus</i>	NCTC <sup>®</sup> 7464 / ATCC <sup>®</sup> 10876	MM86	
<i>Bacillus subtilis</i>	NCTC <sup>®</sup> 10400 / ATCC <sup>®</sup> 6633	MM29	00003
<i>Bacteroides fragilis</i>	NCTC <sup>®</sup> 9343 / ATCC <sup>®</sup> 25285	MM44	
<i>Campylobacter jejuni</i>	NCTC <sup>®</sup> 11351 / ATCC <sup>®</sup> 33560	MM36	
<i>Campylobacter jejuni</i>	NCTC <sup>®</sup> 11322 / ATCC <sup>®</sup> 29428	MM82	00156
<i>Candida albicans</i>	NCPF <sup>®</sup> 3255 / ATCC <sup>®</sup> 2091	MM28	00055
<i>Candida albicans</i>	NCPF <sup>®</sup> 3179 / ATCC <sup>®</sup> 10231	MM42	00054
<i>Citrobacter freundii</i>	NCTC <sup>®</sup> 9750 / ATCC <sup>®</sup> 8090	MM27	
<i>Clostridium perfringens</i>	NCTC <sup>®</sup> 8237 / ATCC <sup>®</sup> 13124	MM45	00007
<i>Clostridium sporogenes</i>	NCTC <sup>®</sup> 532 / ATCC <sup>®</sup> 19404	MM31	00008
<i>Enterobacter aerogenes</i>	NCTC <sup>®</sup> 10006 / ATCC <sup>®</sup> 13048	MM26	00175
<i>Enterobacter cloacae</i>	NCTC <sup>®</sup> 13380 / ATCC <sup>®</sup> 23355	MM01	00082
<i>Enterobacter cloacae</i>	NCTC <sup>®</sup> 13406	MM51	
<i>Enterococcus faecalis</i>	NCTC <sup>®</sup> 775 / ATCC <sup>®</sup> 19433	MM17	00009
<i>Enterococcus faecalis</i>	NCTC <sup>®</sup> 12697 / ATCC <sup>®</sup> 29212	MM18	00087
<i>Enterococcus faecalis</i>	NCTC <sup>®</sup> 13379 / ATCC <sup>®</sup> 51299	MM52	00085
<i>Enterococcus hirae</i>	NCTC <sup>®</sup> 13383 / ATCC <sup>®</sup> 10541	MM35	00011
<i>Escherichia coli</i>	NCTC <sup>®</sup> 12241 / ATCC <sup>®</sup> 25922	MM02	00013
<i>Escherichia coli</i>	NCTC <sup>®</sup> 11954 / ATCC <sup>®</sup> 35218	MM24	
<i>Escherichia coli</i>	NCTC <sup>®</sup> 10418 / ATCC <sup>®</sup> 10536	MM33	
<i>Escherichia coli</i>	NCTC <sup>®</sup> 12923 / ATCC <sup>®</sup> 8739	MM38	00012
<i>Escherichia coli</i>	NCTC <sup>®</sup> 11560	MM63	
<i>Escherichia coli</i>	NCTC <sup>®</sup> 9001 / ATCC <sup>®</sup> 11775	MM75	00090
<i>Escherichia coli</i> CRE	NCTC <sup>®</sup> 13476	MM57	
<i>Escherichia coli</i> (mcr-1)	NCTC <sup>®</sup> 13846	MM34	
<i>Escherichia coli</i> O157 (non-toxigenic)	NCTC <sup>®</sup> 12900 / ATCC <sup>®</sup> 700728	MM93	00014
<i>Haemophilus influenzae</i>	NCTC <sup>®</sup> 8468 / ATCC <sup>®</sup> 9334	MM100	
<i>Haemophilus influenzae</i>	NCTC <sup>®</sup> 12975 / ATCC <sup>®</sup> 49766	MM37	
<i>Haemophilus influenzae</i>	NCTC <sup>®</sup> 12699 / ATCC <sup>®</sup> 49247	MM81	
<i>Haemophilus influenzae</i>	NCTC <sup>®</sup> 11931	MM98	
<i>Klebsiella aerogenes</i>	NCTC <sup>®</sup> 9528	MM88	
<i>Klebsiella pneumoniae</i>	NCTC <sup>®</sup> 9633 / ATCC <sup>®</sup> 13883	MM04	00097
<i>Klebsiella pneumoniae</i>	NCTC <sup>®</sup> 13368 / ATCC <sup>®</sup> 700603	MM83	
<i>Klebsiella pneumoniae</i> CRE	NCTC <sup>®</sup> 13440	MM55	
<i>Klebsiella pneumoniae</i> CRE	NCTC <sup>®</sup> 13443	MM56	
<i>Klebsiella pneumoniae</i> CRE	NCTC <sup>®</sup> 13438	MM58	

## Selectrol Strain Index

Strain Name	Designation	Code	WDCM
<i>Klebsiella pneumoniae</i> CRE	NCTC <sup>®</sup> 13442	MM59	
<i>Lactobacillus brevis</i>	NCTC <sup>®</sup> 13386 / ATCC <sup>®</sup> 8287	MM76	
<i>Legionella pneumophila</i> serogroup 1	NCTC <sup>®</sup> 11192 / ATCC <sup>®</sup> 33152	MM08	00107
<i>Listeria innocua</i>	NCTC <sup>®</sup> 11288 / ATCC <sup>®</sup> 33090	MM92	00017
<i>Listeria monocytogenes</i>	NCTC <sup>®</sup> 7973 / ATCC <sup>®</sup> 35152	MM48	00109
<i>Listeria monocytogenes</i>	NCTC <sup>®</sup> 13372 ATCC <sup>®</sup> 7644	MM77	
<i>Listeria monocytogenes</i>	NCTC <sup>®</sup> 11994	MM87	00019
<i>Neisseria gonorrhoeae</i>	NCTC <sup>®</sup> 8375 / ATCC <sup>®</sup> 19424	MM05	
<i>Neisseria gonorrhoeae</i>	NCTC <sup>®</sup> 12700 / ATCC <sup>®</sup> 49226	MM96	
<i>Proteus mirabilis</i>	NCTC <sup>®</sup> 13376 / ATCC <sup>®</sup> 14153	MM43	
<i>Proteus mirabilis</i>	NCTC <sup>®</sup> 10975	MM68	
<i>Proteus vulgaris</i>	NCTC <sup>®</sup> 4175 / ATCC <sup>®</sup> 13315	MM09	
<i>Pseudomonas aeruginosa</i>	NCTC <sup>®</sup> 12903 / ATCC <sup>®</sup> 27853	MM10	00025
<i>Pseudomonas aeruginosa</i>	NCTC <sup>®</sup> 12924 / ATCC <sup>®</sup> 9027	MM40	00026
<i>Pseudomonas aeruginosa</i>	NCTC <sup>®</sup> 13359 / ATCC <sup>®</sup> 15442	MM41	
<i>Pseudomonas aeruginosa</i>	NCTC <sup>®</sup> 10662 / ATCC <sup>®</sup> 25668	MM65	00114
<i>Rhodococcus equi</i>	NCTC <sup>®</sup> 1621 / ATCC <sup>®</sup> 6939	MM97	00028
<i>Saccharomyces cerevisiae</i>	NCTC <sup>®</sup> 10716/ ATCC <sup>®</sup> 9763	MM50	00058
<i>Saccharomyces cerevisiae</i>	NCPF <sup>®</sup> 3178	MM73	
<i>Salmonella</i> Nottingham	NCTC <sup>®</sup> 7832	MM84	
<i>Salmonella</i> Poona	NCTC <sup>®</sup> 4840	MM89	
<i>Salmonella</i> Typhimurium	NCTC <sup>®</sup> 12023/ ATCC <sup>®</sup> 14028	MM11	00031
<i>Serratia marcescens</i>	NCTC <sup>®</sup> 13382 / ATCC <sup>®</sup> 8100	MM12	
<i>Staphylococcus aureus</i>	NCTC <sup>®</sup> 12981 / ATCC <sup>®</sup> 25923	MM13	00034
<i>Staphylococcus aureus</i>	NCTC <sup>®</sup> 12973 / ATCC <sup>®</sup> 29213	MM14	00131
<i>Staphylococcus aureus</i>	NCTC <sup>®</sup> 7447 / ATCC <sup>®</sup> 6538P	MM30	00033
<i>Staphylococcus aureus</i>	NCTC <sup>®</sup> 10788 / ATCC <sup>®</sup> 6538	MM46	00032
<i>Staphylococcus aureus</i>	NCTC <sup>®</sup> 6571 / ATCC <sup>®</sup> 9144	MM85	00035
<i>Staphylococcus aureus</i> (MRSA)	NCTC <sup>®</sup> 12493	MM64	00212
<i>Staphylococcus aureus</i> (MRSA)	NCTC <sup>®</sup> 13373 / ATCC <sup>®</sup> 43300	MM91	00211
<i>Staphylococcus epidermidis</i>	NCTC <sup>®</sup> 13360 / ATCC <sup>®</sup> 12228	MM15	00036
<i>Streptococcus agalactiae</i>	NCTC <sup>®</sup> 8181 / ATCC <sup>®</sup> 13813	MM16	
<i>Streptococcus pneumoniae</i>	NCTC <sup>®</sup> 12695 / ATCC <sup>®</sup> 6303	MM19	
<i>Streptococcus pneumoniae</i>	NCTC <sup>®</sup> 12977 / ATCC <sup>®</sup> 49619	MM95	
<i>Streptococcus pyogenes</i>	NCTC <sup>®</sup> 12696 / ATCC <sup>®</sup> 19615	MM20	
<i>Vibrio parahaemolyticus</i>	NCTC <sup>®</sup> 10885	MM06	00185
<i>Yersinia enterocolitica</i>	NCTC <sup>®</sup> 12982 / ATCC <sup>®</sup> 9610	MM80	00038

## Selectrol Strains Listed by WDCM Number

WDCM	Strain Name	Designation	Code
00001	<i>Bacillus cereus</i>	NCTC <sup>®</sup> 10320 / ATCC <sup>®</sup> 9634	MM21
00003	<i>Bacillus subtilis</i>	NCTC <sup>®</sup> 10400 / ATCC <sup>®</sup> 6633	MM29
00007	<i>Clostridium perfringens</i>	NCTC <sup>®</sup> 8237 / ATCC <sup>®</sup> 13124	MM45
00008	<i>Clostridium sporogenes</i>	NCTC <sup>®</sup> 532 / ATCC <sup>®</sup> 19404	MM31
00009	<i>Enterococcus faecalis</i>	NCTC <sup>®</sup> 775 / ATCC <sup>®</sup> 19433	MM17
00011	<i>Enterococcus hirae</i>	NCTC <sup>®</sup> 13383 / ATCC <sup>®</sup> 10541	MM35
00012	<i>Escherichia coli</i>	NCTC <sup>®</sup> 12923 / ATCC <sup>®</sup> 8739	MM38
00013	<i>Escherichia coli</i>	NCTC <sup>®</sup> 12241 / ATCC <sup>®</sup> 25922	MM02
00014	<i>Escherichia coli</i> O157 (non-toxigenic)	NCTC <sup>®</sup> 12900 / ATCC <sup>®</sup> 700728	MM93
00017	<i>Listeria innocua</i>	NCTC <sup>®</sup> 11288 / ATCC <sup>®</sup> 33090	MM92
00019	<i>Listeria monocytogenes</i>	NCTC <sup>®</sup> 11994	MM87
00025	<i>Pseudomonas aeruginosa</i>	NCTC <sup>®</sup> 12903 / ATCC <sup>®</sup> 27853	MM10
00026	<i>Pseudomonas aeruginosa</i>	NCTC <sup>®</sup> 12924 / ATCC <sup>®</sup> 9027	MM40
00028	<i>Rhodococcus equi</i>	NCTC <sup>®</sup> 1621 / ATCC <sup>®</sup> 6939	MM97
00031	<i>Salmonella</i> Typhimurium	NCTC <sup>®</sup> 12023 / ATCC <sup>®</sup> 14028	MM11
00032	<i>Staphylococcus aureus</i>	NCTC <sup>®</sup> 10788 / ATCC <sup>®</sup> 6538	MM46
00033	<i>Staphylococcus aureus</i>	NCTC <sup>®</sup> 7447 / ATCC <sup>®</sup> 6538P	MM30
00034	<i>Staphylococcus aureus</i>	NCTC <sup>®</sup> 12981 / ATCC <sup>®</sup> 25923	MM13
00035	<i>Staphylococcus aureus</i>	NCTC <sup>®</sup> 6571 / ATCC <sup>®</sup> 9144	MM85
00036	<i>Staphylococcus epidermidis</i>	NCTC <sup>®</sup> 13360 / ATCC <sup>®</sup> 12228	MM15
00038	<i>Yersinia enterocolitica</i>	NCTC <sup>®</sup> 12982 / ATCC <sup>®</sup> 9610	MM80
00053	<i>Aspergillus brasiliensis</i>	NCPF <sup>®</sup> 2275 / ATCC <sup>®</sup> 16404	MM94
00054	<i>Candida albicans</i>	NCPF <sup>®</sup> 3179 / ATCC <sup>®</sup> 10231	MM42
00055	<i>Candida albicans</i>	NCPF <sup>®</sup> 3255 / ATCC <sup>®</sup> 2091	MM28
00058	<i>Saccharomyces cerevisiae</i>	NCTC <sup>®</sup> 10716 / ATCC <sup>®</sup> 9763	MM50
00082	<i>Enterobacter cloacae</i>	NCTC <sup>®</sup> 13380 / ATCC <sup>®</sup> 23355	MM01
00085	<i>Enterococcus faecalis</i>	NCTC <sup>®</sup> 13379 / ATCC <sup>®</sup> 51299	MM52
00087	<i>Enterococcus faecalis</i>	NCTC <sup>®</sup> 12697 / ATCC <sup>®</sup> 29212	MM18
00090	<i>Escherichia coli</i>	NCTC <sup>®</sup> 9001 / ATCC <sup>®</sup> 11775	MM75
00097	<i>Klebsiella pneumoniae</i>	NCTC <sup>®</sup> 9633 / ATCC <sup>®</sup> 13883	MM04
00107	<i>Legionella pneumophila</i> serogroup 1	NCTC <sup>®</sup> 11192 / ATCC <sup>®</sup> 33152	MM08
00109	<i>Listeria monocytogenes</i>	NCTC <sup>®</sup> 7973 / ATCC <sup>®</sup> 35152	MM48
00114	<i>Pseudomonas aeruginosa</i>	NCTC <sup>®</sup> 10662 / ATCC <sup>®</sup> 25668	MM65
00131	<i>Staphylococcus aureus</i>	NCTC <sup>®</sup> 12973 / ATCC <sup>®</sup> 29213	MM14
00156	<i>Campylobacter jejuni</i>	NCTC <sup>®</sup> 11322 / ATCC <sup>®</sup> 29428	MM82
00175	<i>Enterobacter aerogenes</i>	NCTC <sup>®</sup> 10006 / ATCC <sup>®</sup> 13048	MM26
00185	<i>Vibrio parahaemolyticus</i>	NCTC <sup>®</sup> 10885	MM06
00211	<i>Staphylococcus aureus</i> (MRSA)	NCTC <sup>®</sup> 13373 / ATCC <sup>®</sup> 43300	MM91
00212	<i>Staphylococcus aureus</i> (MRSA)	NCTC <sup>®</sup> 12493	MM64

# Notes

[www.tcsbiosciences.co.uk](http://www.tcsbiosciences.co.uk)

**TCS Biosciences Ltd**

Botolph Claydon, Buckingham, MK18 2LR, United Kingdom

t: +44 (0) 1296 714222, f: +44 (0) 1296 714806, e: [sales@tcsgroup.co.uk](mailto:sales@tcsgroup.co.uk)

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