


ECOLOGICAL THERMOMETERS - MERCURY FREE

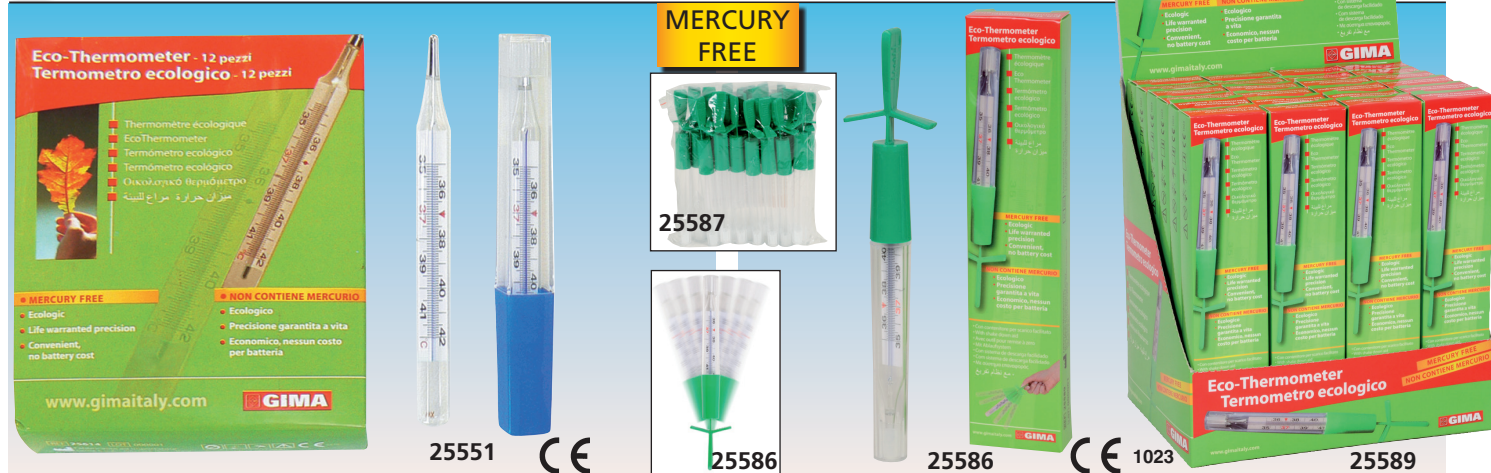

- 25551 ECOLOGICAL OVAL THERMOMETER - box of 12
 - 25587 SHAKER for ecological thermometers - box of 24
- Mercury free clinical thermometer with gallium alloy instead of mercury. Safe, precise and functional. Glass structure. Multilingual box and instructions: GB, FR, IT, ES, PT, DE, PL, RO, SK, GR, Arabic.

- 25586 ECOLOGICAL OVAL THERMOMETER with shaker - box of 1
 - 25589 ECOLOGICAL OVAL THERMOMETER with shaker - display box of 24 pcs.
- Same as 25551 but supplied with a shaker to help shake down gallium after use.

WIRELESS THERMOMETER BLUETOOTH 4.0

- 25548 TEMP SITTER - baby wireless thermometer
- Accurate and innovative thin thermometer. Transmission to Temp Sitter App (iOS, Android). To be placed under baby armpit. Update frequency every 3 seconds. High temp alarm can be set in the App. Waterproof and material of safe quality. Supplied with 10 adhesives to fix thermometer to baby armpit, manual (GB, FR, IT, ES).



25548

TECHNICAL SPECIFICATIONS

Battery life: 3,000 hours, not replaceable
Measurement range: 25-45 °C (± 0.05 between 35-38.5 °C)
Size: 51x31.8x8.5 mm - Weight: 9.6 g - Switchable °C/°F

CE 0123

Download on the App Store
GET IT ON Google play

OMRON DIGITAL THERMOMETERS

- 25603 ECO TEMP BASIC THERMOMETER - MC 246-E
- 25604 FLEX TEMP SMART THERMOMETER - MC 343F-E




Model	TECHNICAL SPECIFICATIONS	
	25603 (Eco Basic)	25604 (Flex Smart)
Measurement Type	Oral, Armpit, Rectal	Oral, Armpit, Rectal
Flexible tip		•
Beeper (signals when a measurement is complete)	•	•
Memory	last reading	last reading
°C/°F selection	•	•
Large, 3 digit display	•	•
Water resistant	•	•
Automatic shut off	•	•

LIQUID CRYSTAL FOREHEAD THERMOMETERS

- 25599 FEVERTEST - bag of 10 blister
- 25596 FEVERTEST BULK - box of 100 - hospital use
- 25600 FEVERCONTROL - bag of 10 blister

Reusable, non toxic, mercury free, safe, unbreakable, reliable, long lasting, easy to read, liquid crystal thermometers. They indicate temperature in 15 sec with graduation 0.5°C. Medical Devices.

MADE IN ITALY

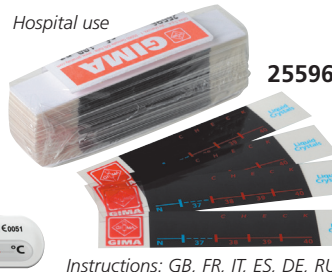

Instructions: GB, FR, IT, ES, DE, RU.

48h FEVERMONITOR

- 25601 48h FEVER MONITOR - bag of 10 blister of 3 pcs.

Innovative "continual" fever indicator for 48h monitoring, for babies. Disposable, non toxic, mercury free, safe, unbreakable, reliable, easy to read, liquid crystal thermometers. It indicates temperature in 30 seconds. Medical Device.



MADE IN ITALY

Blister in GB, IT

CE 0051

Instruction For Use

a. Mercury free clinical thermometer is a kind of measuring device for patient's temperature, and it can be used to measure the temperature of patient's mouth or armpit or anus.

b. The product is in compliance with EN12470-1:2000+A1:2009, Clinical thermometers - Part 1: Metallic liquid-in glass thermometers with maximum device.

1. The conditions of use: measuring temperature range should be between 35 °C to 42 °C, and the environment temperature should be 20 °C to 30 °C; if the temperature for storage is less than 20 °C, please put the thermometers in the temperature 20 °C to 30 °C for at least 1-2 hours or 30 °C to 35 °C for ten minutes before use or test.

2. The temperature must be less than 42 °C when it is stored and transported;

3. Storage of the thermometer: after use, swing the metallic liquid down.

This procedure will require a greater effort in comparison with a common mercury thermometer.

Then clean it and put it plastic case to avoid any damage.

c. Operation procedures

1. The liquid column shall fall below the lowest numbered scale mark before use, then clean and sterilize the clinical thermometer in medical alcohol or put on a sterilized plastic tube.

2. Be careful to put it below the tongue of patients for oral clinical thermometer and tell them to close their mouth tightly.

3. A clinical thermometer for rectal use is applied some lubricator before inserting it in the anus.

4. Be careful to put it below the armpit of patients if measuring armpit position and tell them to tight it tightly.

Please read the number of the clinical thermometer about 5 minutes after above measuring steps.

d. Warning

1. It is not suitable to use oral clinical thermometers for not more than ten years' children;

2. Warn patients against chewing, speaking and playing when the clinical thermometers are below the tongue of patients;

3. Need to let sick children lie quietly and help them to hold another end of the clinical thermometers for rectal use in case it is inserted too deep, or it slips from the anus.












4. Avoid having a bath, eating or smoking half an hour before taking body temperature;

5. It is not designed for special applications (e.g. thermometers for premature babies, ovulation thermometers);

6. Inspect the clinical thermometer once a year.


e. Disposal of the Clinical thermometers

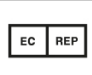The metallic liquid is an alloy consisting of gallium, indium and tin. this clinical thermometer should be handled with care. In particular, care should be taken to avoid dropping the instrument or treating it in any way that could result in damage. When the clinical thermometers are broken, cleaning split metallic liquid by a damp cloth with soap suds.

	Product code		Keep away from sunlight		Medical Device complies with Directive 93/42/EEC
	Lot number		Keep in a cool, dry place		Caution: read instructions (warnings) carefully
	Manufacturer		Consult instructions for use		Authorized representative in the European community
	Date of manufacture		Mercury free		



Via Marconi, 1 - 20060 Gessate (MI) Italy
gima@gimaitaly.com - export@gimaitaly.com
www.gimaitaly.com

 25586 / CR.W00

M 25550 - EN - Rev. 2.02.20

 Wuxi Medical Instrument Factory Co., Ltd.
No. 43 Xixin Road, ZhangJing, Xibei Town,
Wuxi city, Jiangsu, 214194, China
Made in China.



 Lotus NL B. V.
Koningin Julianaplein 10, 1e Verd, 2595 AA,
The Hauge, Netherlands.