



Laboratory Equipment Manufacturer
www.mrclab.com



Operation Manual for BIOLOGICAL MICROSCOPE BAM200 Series



This operation manual is for the biological microscope of BAM200 series. A thorough and careful reading of the manual is advised for safe and easy use of the microscope.

3, Hagavish st. Israel 58817 Tel: 972 3 5595252, Fax: 972 3 5594529 mrc@mrclab.com

MRC.VER.01-11.11

Contents

Correct assembly and adjustments are critical for the microscope to exhibit its full performance. If you are going to assemble the microscope yourself, please read Chapter2, "ASSEMBLY" (pages6-8) carefully. For the modules provided with instruction manuals, also read the assembly procedures in their instruction manuals.

IMPORTAT

Be sure to read this section for safe use of the equipment...1-3

1 SAFETY SYMBOLS ON THE MICROSCOPE.....	1
2 Caution Symbols in this Manual.....	1
3 Safety Precautions	1
4 Moving the Microscope.....	2
5 Working environment	3
6 Electric Power Specifications.....	3
7 Maintenance and Storage.....	3
8 Unpacking the Microscope.....	3

1 MODULE & CONTROLS NOMENCLATUE..... 4-5

2 ASSEMBLY..... 6-8

2.1 Assembly Diagram.....	6
2.2 Detailed Assembly Procedures	6
2.2.1 Attaching/Replacing the Bulb.....	6
2.2.2 Attaching the Observation Tube.....	7
2.2.3 Attaching the Eyepieces.....	7
2.2.4 Attaching the Objectives	7
2.2.5 Attaching the Condenser	8
2.2.6 Mounting the Filter.....	8

3 LIGHT BRIGHTFIELD OBSERVATION..... 9-13

3.1 Light Brightfield Observation Procedure.....	9
3.2 Using the Controls.....	10
3.2.1 Power Supply Unit.....	10
3.2.2 Focusing Block.....	10
a) Rotation Direction of the Coarse/Fine Adjustment Knobs.....	10
b) Adjusting the Coarse Adjustment Knob Tension.....	10
3.2.3 Mechanical Stage.....	10
a) Placing the Specimen on the stage.....	10

b) Moving the Specimen.....	10
c) Raise and lower the stage.....	10
3.2.4 Observation Tube.....	10
a) Adjusting the Interpupillary Distance.....	10
b) Adjusting the Diopter.....	11
c) Using the Eye Shades.....	11
3.2.5 Condenser.....	11
a) Centering the Condenser.....	11
b) Using the Aperture Iris Diaphragm.....	12
c) Condenser height adjustment.....	12
3.2.6 Objectives.....	12
a) Engage the Required Objectives.....	12
b) Using Immersion Objectives.....	12
3.2.7 Filters.....	13

4 OTHER OBSERVATION METHODS..... 14-15

4.1 Phase Contrast Observation.....	14
4.2 Simplified Polarized Light Observation.....	14
4.3 Dark Field Light Observation.....	15
4.4 CCD & Computer Observation.....	15
4.5 Digital Camera Observation.....	15
4.6 Camera.....	15

5 TROUBLESHOOTING GUIDE..... 16-18

6 STANDARD OUTFITS & SPECIFICATIONS..... 19-22

IMPORTANT

1 SAFETY SYMBOLS ON THE MICROSCOPE

The following symbols are found on the microscope. Study the meaning of the symbols and always use the equipment in the safest possible manner.

Symbol	Explanation
	High temperature warning! Indicates that the surface becomes hot, and should not be touched with bare hands. This symbol is on the bottom cover of microscope.
	Beside the fuse and power connector. Do not touch! Electric shock may happen when the power cord is connected.
I	Indicates that the main switch is on.
O	Indicates that the main switch is OFF.

2 Caution Symbols in this Manual

If the microscope is used in a manner not specified by this manual, the safety of the user may be imperiled. In addition, the equipment may also be damaged. Always use the equipment as outlined in this instruction manual. The following symbols are used to set off text in this instruction manual.

Symbol	Explanation
	Indicates that failure to follow the instructions in the warning could result in bodily harm to the user and /or damage to equipment (including objects in the vicinity of the equipment).
★	Indicates that failure to follow the instructions could result in damage to equipment.
©	Indicates commentary (for ease of operation and proper maintenance).

3 Safety Precautions

3.1 After the equipment has been used in an observation of a specimen that is accompanied with a potential of infection, clean the parts coming in contact with the specimen to prevent infection.

- Moving this product is accompanied with the risk of dropping the specimen. Be sure to remove the specimen before moving this product.
- In case the specimen is damaged by erroneous, promptly take the infection prevention measures.

3.2 The microscope is provided with a simplified waterproof mechanism. Therefore, if culture liquid or water is split on the stage, revolving nosepiece or microscope frame, damage to the equipment or an electrical shock may result. Immediately wipe the liquid or water off if it is spilt on them.

3.3 The microscope is not covered by warranty in terms of laser safety. The user should assume liabilities for any consequence of user modification,

3.4 The surfaces of the lamp housing will become extremely hot during long-time operation. Be sure to keep the flammable stuffs such as paper, alcohol, oil away from the lamp house to avoid fire.

3.5 When using the microscope, route the power cord away from the lamp housing. Should the power cord come in contact with the hot lamp housing, the power cord could melt and cause electric shock.

3.6 To avoid potential shock hazards and burns when replacing the light bulb, set the main Switch to “O” (OFF) then disconnect the power cord from the wall outlet in advance. Whenever you replace the bulb during use or right after use, allow the lamp housing and bulb to cool before touching.

3.7 Electric shock warning:

Remove of the bottom cover of the microscope makes the dangerous electric parts inside exposed. Any contact with these parts may cause shock or death. In event of maintenance, please apply to qualified professionals for help.

3.8 The G4 bulb socket is designed specially for 6V/20W halogen bulb. Damage will occur if bulb of different description is replaced.

3.9 Always be sure the power cord provided by UOP. If the proper power cord is not used, product safety performance cannot be warranted.

3.10 To avoid potential shock hazards when replacing the fuse, set the main switch to “O”(OFF) then disconnect the power cord from the wall outlet in advance.

3.11 Always ensure that the grounding terminal of the microscope and that of the wall outlet are properly connected. If the equipment is not grounded, UOP can no longer warrant the electrical safety performance of the equipment.

3.12 Never insert metallic objects into the air vents(A^⑨ in Fig.2) of the microscope frame as this could result in electrical shock, personal injury and equipment damage.

3.13 A microscope is a precision instrument. Handle it with care and avoid subjecting it to sudden or severe impact.

4 Moving the Microscope

4.1 When moving the microscope, remove the observation tube, then carefully carry the microscope frame by the base (front edge)(^② in Fig.1) and the

grasping part(^① in Fig.1) in the upper rear.

4.2 Also be sure to remove the specimen since it may fall.

4.3 When moving the microscope for a long distance, it is also recommended to disconnect all cables from the equipment.

4.4 When transporting it, also engage the adhesive tape lock mechanisms and package it sufficiently.

4.5 Also be careful against slipping of hands during carrying.

★ Damage to the microscope will occur if you grasp it by other parts including the stage, coarse/fine adjustment knobs, the nosepieces, etc.

5 Working environment

5.1 Indoor use.

5.2 Ambient Temperature: 5°C to 40°C (41°F to 104°F)

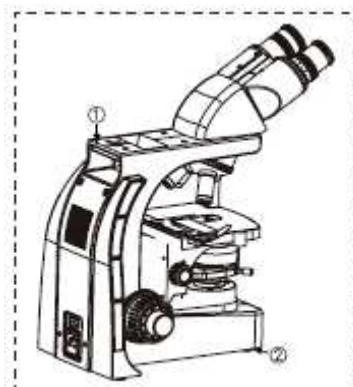


Fig.1

5.3 Maximum relative humidity: 80% for temperatures up to 31°C (93°F), 60% at 37°C (99°F), to 50% relative humidity at 40°C (104°F).

5.4 Supply voltage fluctuations: $\pm 10\%$

5.5 Pollution degree: 2 (In accordance with IEC60664)

6 Electric Power Specifications

Input: 100-240V ~ 0.5A, 47-63Hz

Output: 6V 3.4A

Fuse: 3.15A, 250V, F $\phi 5 \times 20$ mm

Halogen Bulb Socket: G4

7 Maintenance and Storage

7.1 Clean all glass components by wiping gently with gauze. To remove fingerprints of oil smudges, wipe with gauze slightly moistened with a mixture of ether (70%) and alcohol (30%). Since solvents such as ether and alcohol are highly flammable, they must be handled carefully. Be sure to keep these chemicals always from open flames, or potential sources of electrical sparks—for example, electrical equipment that is being switched on or off. Also remember to always use these chemicals only in a well-ventilated room.

7.2 Be sure to clean the oil immersion objective after use. Leaving immersion oil on it will degrade its performance.

7.3 Do not attempt to use organic solvents to clean the non-optical components of microscope. To clean them, use a lint-free, soft cloth slightly moistened with a diluted neutral detergent.

7.4 Never attempt to disassemble any part of the microscope.

7.5 When not using the microscope, make sure to set the main switch to "O" (OFF), confirm that the lamp housing is cool enough and cover the microscope with the provided dust cover.

7.6 Do not use the microscope where it is subjected to direct sunlight, high temperature and humidity, dust or vibrations.

8 UNPACKING THE MICROSCOPE

8.1 Please check all the components according to the packing list in the package as you unpacking the microscope. Contact us or our distributor as soon as possible if any component is missed in the package.

8.2 Before transporting the microscope, we have fixed the flexible parts of the stage with pieces of adhesive tape, please remove the adhesive tapes before use.

1. MODULES & CONTROLS NOMENCLATURE

◎ The modules shown below are only the representative modules. As there are other

modules which can be combined with the microscope but are not shown below. Please also refer to the latest UOP catalogues or your dealer.

◉ If you have not yet assembled the microscope, read Chapter 2, “Assembly” (pages6-8)

• The illustration shows the system composed of modules and controls enclosed in

□.

RIGHT SIDE VIEW OF MICROSCOPE

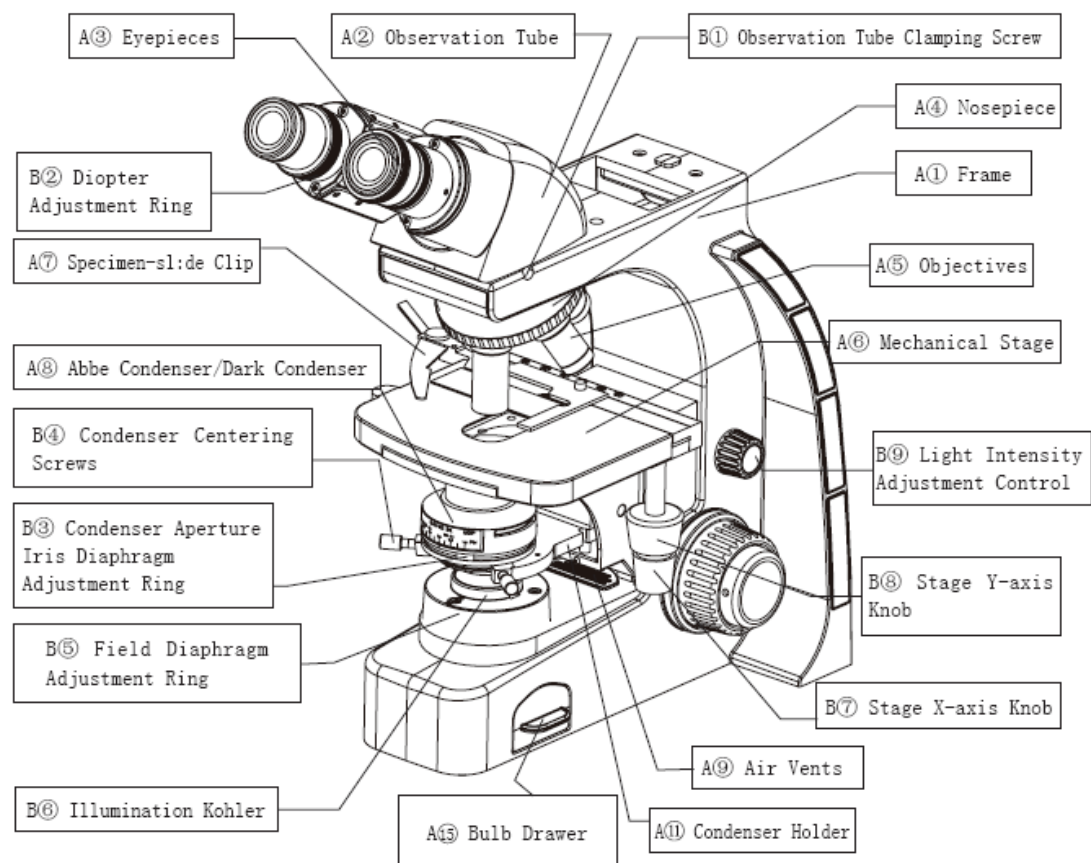


Fig.2

LEFT SIDE VIEW OF MICROSCOPE

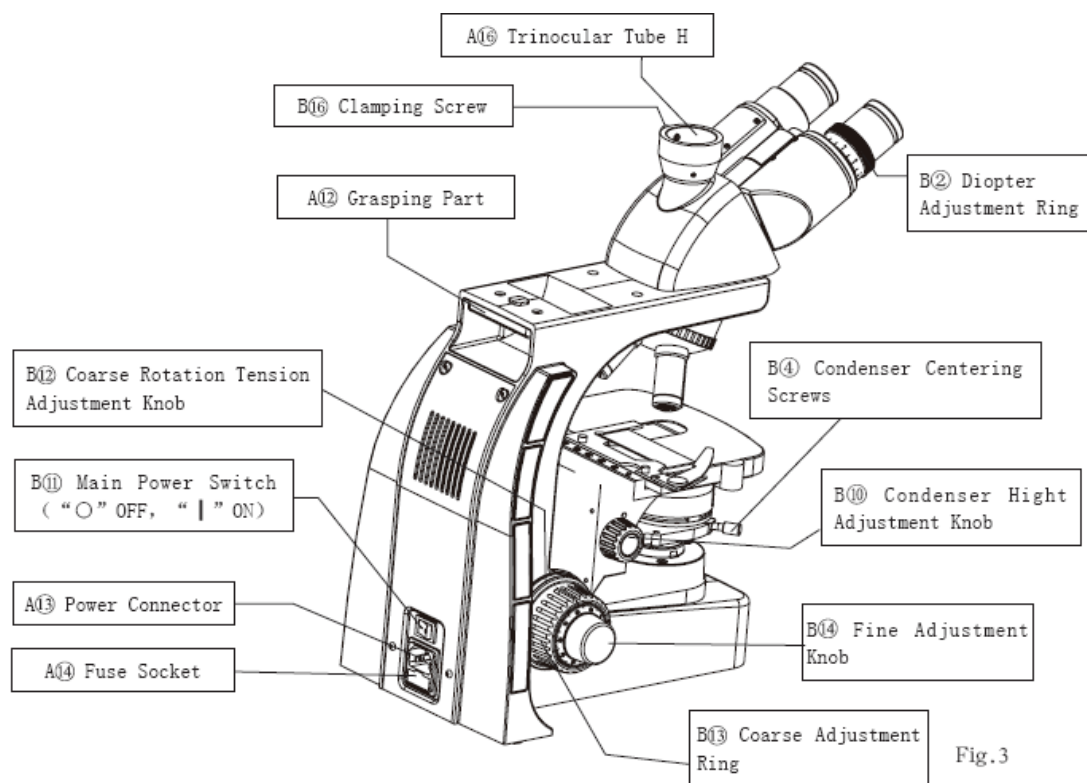


Fig.3

A1 Frame	A2 Observation Tube	A3 Eyepieces
A4 Nosepiece	A5 Objectives	A6 Mechanical Stage
A7 Specimen-slide Clip	A8 Abbe Condenser/Dark Condenser	A9 Air Vents
A10 Filter Holders(Figs.6&7)	A11 Condenser Holder	A12 Grasping Part
A13 Power Connector	A14 Fuse Socket	A15 Bulb Drawer(Fig.5)
A16 Trinocular Tube H		
B1 Observation Tube Clamping Screw	B2 Diopter Adjustment Ring	
B3 Condenser Aperture Iris Diaphragm Adjustment Ring	B4 Condenser Centering Screws(two)	
B5 Field Diaphragm Adjustment Ring	B6 Illumination Kohler	
B7 Stage X-axis Knob	B8 Stage Y-axis Knob	
B9 Light Intensity Adjustment Control	B10 Condenser Height Adjustment Knob	
B11 Main Power Switch ("O" OFF, "I" ON)	B12 Coarse Rotation Tension Adjustment Ring	
B13 Coarse Adjustment Knob	B14 Fine Adjustment Knob	
B15 Grip of the Bulb Drawer	B16 Clamping Screw	

2.ASSEMBLY

2.1 Assembly Diagram

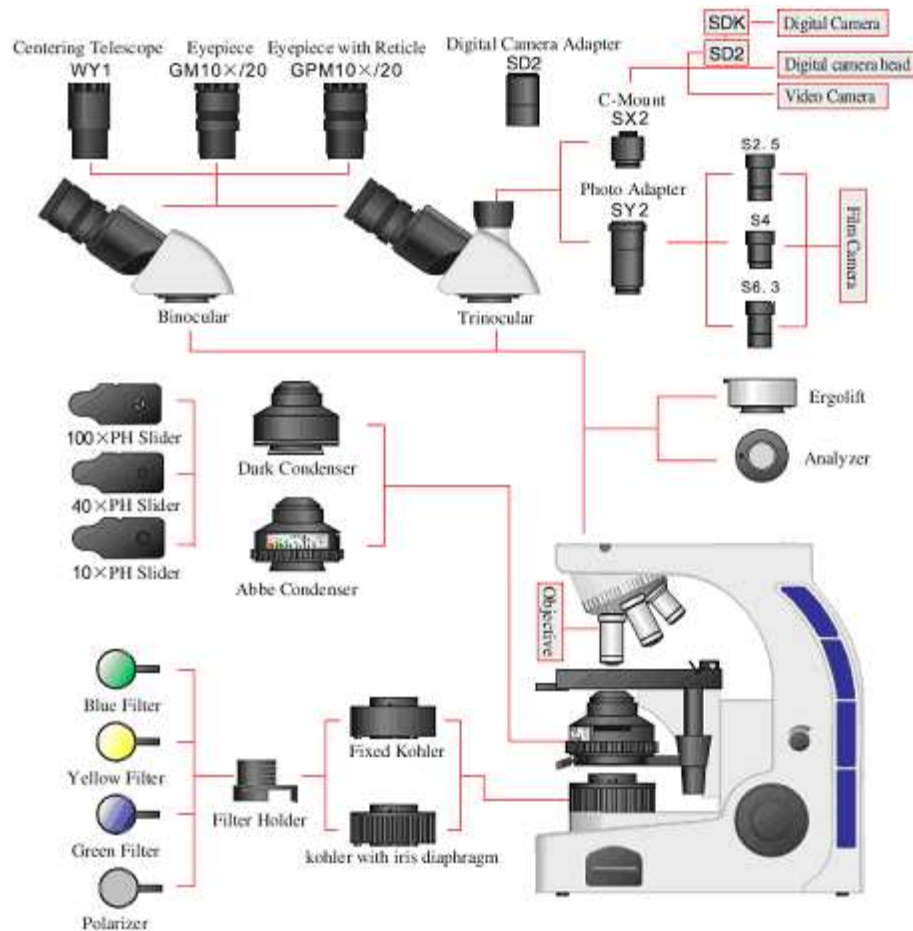


Fig.4

2.2 Detailed Assembly Procedures

2.2.1 Attaching/Replacing the Halogen Bulb(② in Fig.5)

Caution for Bulb Replacement During or Right After Use.

To avoid potential shock hazards and burns when replacing the light bulb, set the main switch(B in Fig.3) to "O" (OFF) . Then disconnect the power cord from the wall outlet in advance. Whenever you replace the bulb during use or right after use, allow the lamp housing and bulb to cool before touching.

To prevent reduced bulb life or cracking, do not touch the bulb with bare hands. If finger prints are accidentally left on the bulb, wipe the bulb with a soft cloth.

a) Grip the Grip(B₁₅ in Fig.5) of the bulb drawer and pull out the bulb drawer (on the right side of the microscope base, A₁₅ in Fig.5).

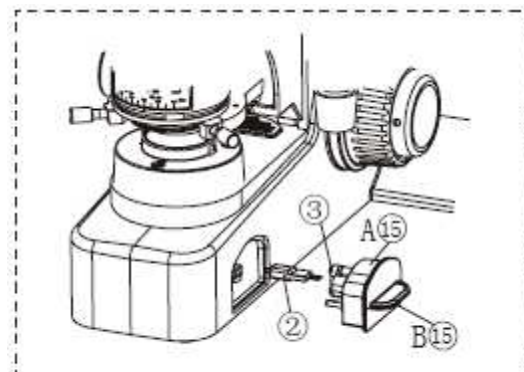


Fig. 5

b) Pull out the burnt-out bulb from the bulb socket if you want to replace it.

d) Hold the new bulb with gloves or a piece of gauze, insert the bulb pins and fully into the pin sections on(③ in Fig.5) the bulb socket.

e) Push in the bulb drawer and close the lamp house.

⊙ After working for above 10 hours, cut off the instrument about 30 minutes to prolong the bulb life as possible.



2.2.2 Attaching the Observation Tube(Figs.2&4)

a) Using the Allen screwdriver, loosen the observation tube clamping screw (B① in Fig.2)on the observation tube mount.

b) Attach the circular dovetail mount of the observation tube into the observation tube mount , placing the observation tube so that the interpupillary distance scale numbers are seen right side up. Then tighten the clamping screw to clamp the observation tube.

2.2.3 Attaching the Eyepieces(Figs.2&4)

a) Remove the eyepieces's dust caps.

b) Insert the eyepiece into the eyepiece sleeve.

2.2.4 Attaching the Objectives(Fig.2)

Be sure no specimen is on the stage before you attaching the objectives to prevent potential damage to the specimen-slide.

a) Lower the stage to the farthest position.

b) Screw the objectives into the nosepiece(A④ in Fig.2) in the order from lower-power to higher-power.

c) Please cover the empty positions in the nosepiece with the dust cap to protect the optical parts from the dust.

2.2.5 Attaching the Condenser(Figs.2&4)

a) Raise the stage to the farthest position by rotating the coarse adjustment knob (B13 in Fig.3) and lower the condenser holder to the farthest position by rotating the condenser height adjustment knob(B⑩ in Fig.3)

b) Loosen the condenser centering screws.(B④ in Fig.2)

c) Fit the required condenser into the mount dovetail on the condenser holder(A11 in Fig.2), and push in the condenser until its positioning pin fits into the positioning

groove on the mount dovetail of the condenser holder.

d) Tighten the condenser centering screws.

e) Centering the condenser.(see page11, 3.2.5.a)

2.2.6 Mounting the Filter(Figs. 4、 6&7)

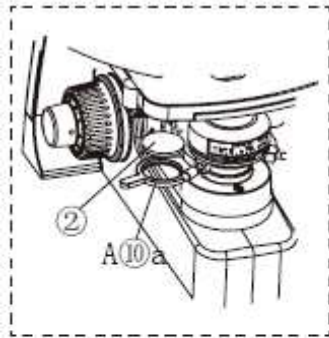


Fig. 6

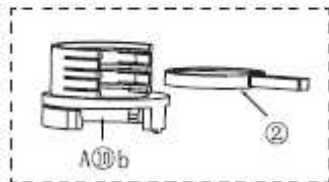


Fig. 7

a) Movable filter holder (attached to the bottom of condenser holder) and filter holder (attached to the illumination kohler) are both available for BAM200 series microscope.

b) Mount the filter to the movable filter holder(Fig.6)

b.1 Stir the filter holder (A10a in Fig.6)out to the left from the bottom of the condenser holder(A11 in Fig.2)

b.2 Insert the filter(2 in Fig.6) down into the filter holder

without tilting. Move the holder(A10a in Fig.6) back to the position right down the condenser into the light path.

c) Mount the filter to the filter holder.(Fig.7)


c.1 Push the filter holder(A10b in Fig.7)downwards and mount it to the illumination Kohler. (B6 in Fig.2)

c.2 Insert the required filter (2 in Fig.7) into the holder. More than one filter can be stacked in this filter holder.

3.LIGHT BRIGHTFIELD OBSERVATION PROCEDURE

3.1 LIGHT BRIGHTFIELD OBSERVATION PROCEDURE

© The following flow shows the operating procedure for the transmitted light brightfield observation which is the basic observation method of this microscope. The operating procedures for phase contrast observation, dark field and Polarization observation will be described separately in Chapter5, "OTHER OBSERVATION METHODS"(Pages 14-15)

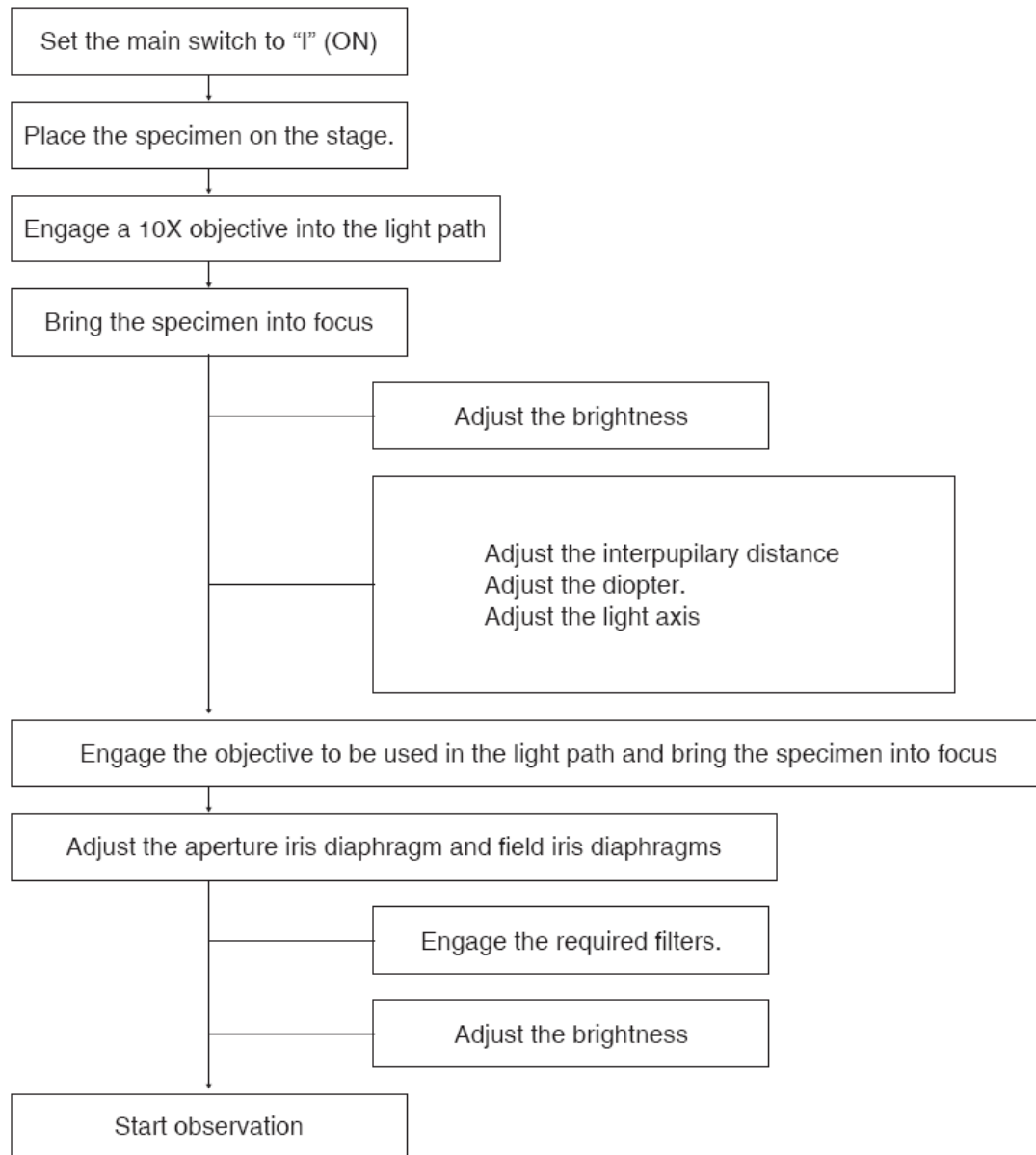

3.2 USING THE CONTROLS

3.2.1 Turning Power On, Adjusting the Brightness(Figs.2&3)

- Make sure that the light intensity control (B[⊙] in Fig.2) is in the MIN (minimum intensity) position and set the main switch(B Fig.3) to "I" (ON).
- Rotate the Light intensity control toward MAX (maximum intensity) to increase the intensity and the illumination brightness.

3.2.2 Focusing Block(Figs.3&8)

- Rotation Direction of the Coarse/Fine Adjustment Knobs

⊙ Rotating the coarse or fine focus adjustment knob(B ,B in Fig.8) toward the front lowers the stage and toward the rear raises the stage.

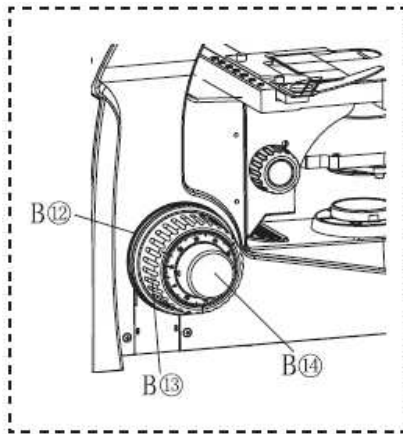


Fig. 8

★ Never rotating the right and left coarse adjustment knobs reversely at the same time, this will damage the focusing block.

b) Adjusting the Coarse Adjustment Knob Tension(Fig.8)

★ Always use the rotation tension adjustment ring(B in Fig.8) to control the rotation tension of the coarse adjustment knob. The tension of the coarse adjustment knob has been pre-adjusted to optimum tension, but this can be changed as required. Turn the rotation tension adjustment ring toward the rear to increase or toward the front to decrease the knob's tension. If the stage lowers by its own weight or the focusing obtained with the fine adjustment knob is lost soon, the tension is set too low. In this case, turn the rotation tension adjustment ring toward the rear to increase the tension.

to increase or toward the front to decrease the knob's tension. If the stage lowers by its own weight or the focusing obtained with the fine adjustment knob is lost soon, the tension is set too low. In this case, turn the rotation tension adjustment ring toward the rear to increase the tension.

3.2.3 Mechanical Stage (Figs.2, 9&10)

a) Placing the Specimen on the stage. Place the specimen on the center of the stage. Attach the stage clips and clamp the specimen

b) Moving the Specimen To move the specimen to a desired position, rotate the X-axis knob and Y-axis knob to move the stage. The travel area is 76mm(X-axis) x 52mm(Y-axis).

c) Raise and lower the stage. (Ref: page 10, 3.2.2.a)



Fig. 9

3.2.4 Observation Tube(Figs.2, 11&12)

a) Adjusting the interpupillary Distance(Fig.11) While looking through the eyepieces, adjust the binocular vision until the left and right fields of view coincide completely.

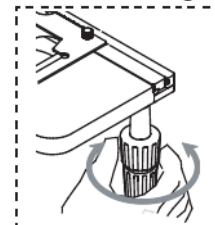


Fig. 10

◎ The interpupillary Distance of BAM200 is 52-74mm.

b) Adjusting the Diopter. (Fig.12)

◎ The diopter adjustment accuracy can be improved by using an objective with as high power as possible.

b.1 While looking through the right eyepiece, rotate the coarse/fine adjustment knobs to bring the specimen into focus.

b.2 Look through the left eyepiece and rotate only the diopter adjustment ring(B② in Fig.2) on the left eyepiece sleeve to bring the specimen into focus.

★ When rotating the diopter adjustment ring of the left eyepiece, hold the lower part of the left eyepiece with the other hand.

c) Using the eye Shades

c.1 When Wearing eyeglasses Use with the eyeshades in the normal, folded-down

position. This will prevent the eyeglasses from being scratched.

c.2 When Not Wearing Eyeglasses Extend the folded eye shades in the upward direction to prevent extraneous light from entering between the eyepieces and eyes.

3.2.5 Condenser (Numerical aperture (N.A.) 1.25 oil)

a) Centering the Condenser(Figs.2, 13,14&15)

⊙ Your microscope illumination Kohler is with field diaphragm or without field diaphragm as you required. The center of condenser had bee accurately centered.

a.1 If you have to re-center the condenser when the microscope illumination Kohler is with field diaphragm (Fig.13)

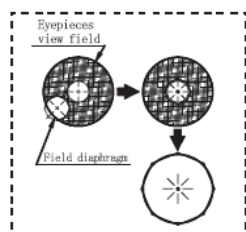


Fig. 13

① Slide the field diaphragm adjusement ring(B⑤ in Fig.2) to the left to open the diaphragm.

② Slide the field iris diaphragm ring to the fully open position.

③ Engage the 10X objective and bring the specimen into focus.

④ Using the field iris diaphragm ring, stop down the field iris diaphragm until its image is just inside the field of view.

⑤ Rotate the condenser height adjustment knob to bring the field iris diaphragm image into focus.

⑥ While gradually opening the field iris diaphragm, rotate the condenser centering

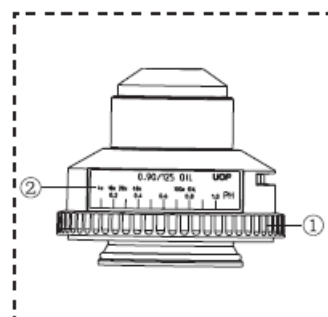


Fig. 14

Screws(B④ in Fig.2) on the condenser holder to adjust so that the field iris diaphragm image is centered in the field of view of the eyepieces.

⑦ To check centration, open the field iris diaphragm until its image inscribed the field of view. Now the condenser is centered.

⊙ In actual observation, open the field iris diaphragm until its image circumscribed the field of view.

a.2 If you have to re-center the condenser when the

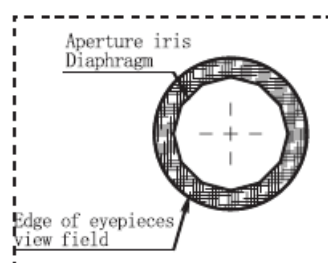


Fig. 15

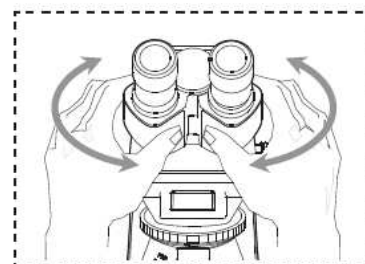


Fig. 11

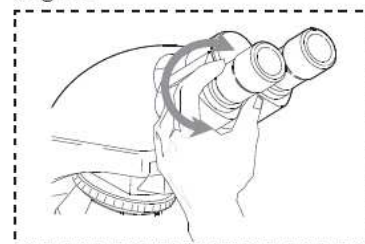


Fig. 12

microscope illumination Kohler is without field diaphragm,

① Stop down the condenser aperture iris diaphragm to correspond with a 10x objective.

② Engage a 10x objective into the light path and bring the specimen into focus.

③ Rotate the condenser centering screws on the condenser holder to adjust so that the specimen image is flat, evenly bright, and full in the field of view of the 10x eyepiece.

④ Now the condenser is centered.

b) Using the Aperture Iris Diaphragm(Fig.14&15)

⊙ In general, the potential resolving power of an objective is fully utilized if the diaphragm is stopped down to correspond with the numerical aperture (N.A.) of the objective.

⊙ Depending on the specimen, image contrast of focal depth in observation or photomicrography may be improved by keeping the aperture iris diaphragm stopped down a little. In general, a good image is obtained if the diaphragm is stopped down to between 70% and 80% of the N.A. of the objective. Stop further down for less contrasty specimens.

⊙ To check the position of the perimeter of the aperture iris diaphragm, remove the eyepieces and look into the eyepiece sleeves to view the aperture iris diaphragm image and the objective's exit pupil. Slide the field iris diaphragm ring(① in Fig.14) to the right to open the diaphragm, and slide the aperture iris diaphragm ring to the left to shut the diaphragm.

⊙ There is color coded mark(② in Fig.14) on the condenser. The magnification of the objectives, the PH objective, the N.A. of the condenser are clearly coded with different color in the mark. For different objective observation, the white line on the iris stop at the corresponding position can bring you best observation.

c) Condenser Height Adjustment: Rotate the condenser height adjustment knob(B⑩ in Fig.3) toward the front to lower the condenser, toward the rear to raise the condenser. Adjust the height of the condenser and to bring the field diaphragm image into focus.

3.2.6 Objectives

a) Rotate the revolving nosepiece by its rim(① in Fig.16) to engage the required objective into the light path and

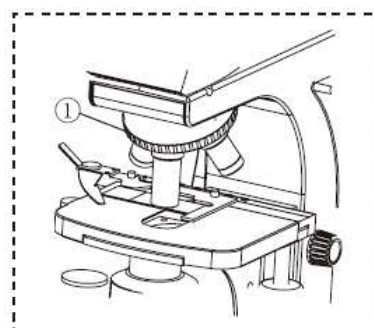


Fig. 16

stop at the “click” position.

b). Do not rotate the nosepiece by pushing the objectives(A⑤ in Fig.2), this will damage the positioning accuracy of the objectives.

c) Using 100X Immersion objectives (Fig.17)

c.1 Using a 40X objective, bring the specimen into focus.

c.2 Rotate the revolving nosepiece counterclockwise, stop down in the position where the specimen is just between the 40X objective and 100X objective.

c.3 Apply a drop of the provided immersion oil to the specimen cover, then rotate the revolving nosepiece to engage the oil immersion objective into the light path.

★ Please remove the oil immediately after use.

★ If the oil contains air bubbles, the image will be degraded. Make sure the oil is free of air bubbles. To remove air bubbles, slightly rock the revolving nosepiece manually to engage and disengage the oil immersion objective once or twice. After use, wipe away the

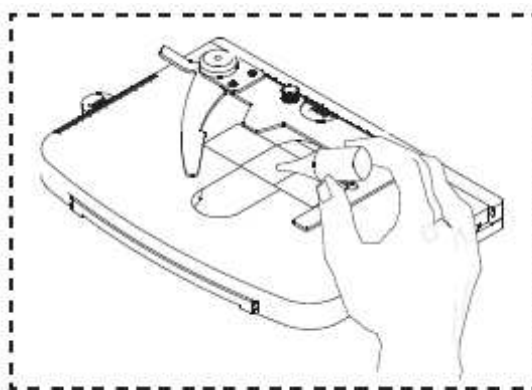


Fig. 17

immersion oil at the objective front lens and the specimen with a piece of gauze lightly moistened with a mixture of ether(70%) and alcohol(30%) Caution on using the immersion oil:

If immersion oil comes into contact with your eye or skin, immediately take the following action.

For eye: Rinse with clean water(for more than 15 minutes),

For skin: Wash with soap and water.

If the appearance of the eye or skin changes or pain continues, immediately consult your doctor.

3.2.7 Filters

a) Filters of different colors are required according to different tinct specimens.

b) If there is no needs of the filter or requiring to enhance the light intensity in the view field while the plan achromatic objective 100×(oil) is being employed, stir the movable filter holder towards the left side out from the light path.

4.OTHER OBSERVATION METHODS

4.1 Phase Contrast Observation(Fig.18)

4.1.1 Phase Contrast Optical Elements and Applicable Objectives

Optical Element	Applicable Objectives
10X Phase slide	10X/0.30 PH
40X Phase slide	40X/0.70 PH

100X Phase slide

100X/1.25(oil) PH

4.1.2 Phase Contrast Observation

a) Conduct the preparation according to the brightfield observation procedure.

b) Slide the aperture iris diaphragm ring(B③ in Fig.18) toward far right to the PH position.

c) Rotate the revolving nosepiece and engage the applied PH objective into the light path.

d) Slowly slide the corresponding Phase slide(① in

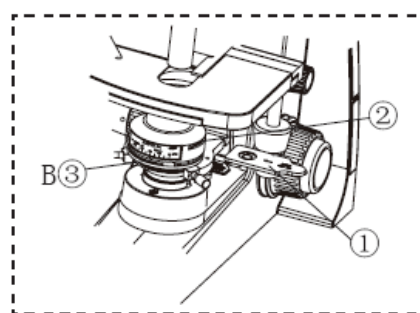


Fig. 18

Fig.18) into the socket(② in Fig.18) on the right side of the condenser (be sure the side of the slide with marks facing upside).

e) Then you are ready to start your phase contrast observation.

4.2 Simplified Polarized Light Observation(Fig.19,20&24)

◎ An analyzer and polarizer are required for simplified polarized light observation.

4.2.1 Attaching the Analyzer and Polarizer.

a) Using the Allen screwdriver, loosen the observation tube clamping screws(B① in Fig.19) on the observation

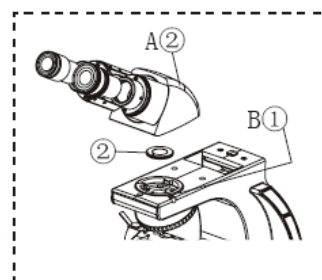


Fig. 19

tube mount and remove the observation tube(A② in Fig.19).

b) Hold the polarizer(② in Fig.19), insert the polarizer positioning pin fully into the pin section of the observation tube mount. Be sure the polarizer is tightly and evenly placed in the observation tube mount.

c) Reassembly the observation tube. (Ref:page 7,2.2.2)

d) Mount the filter holder(① in Fig.20) to the

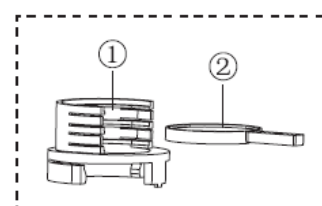


Fig. 20

illumination Kohler.

e) Insert the analyzer(② in Fig.20) into the filter holder with the symbol side upward.

4.2.2 Polarized light observation:

a) Before your polarized light observation, do your preparation as the bright light observation.

b) Rotate the analyzer and stop in the position when your eyepieces view field is of the darkest.

c) Place the required specimen on the stage and start your polarize light observation.

4.3 Dark-Field Light observation

4.3.1 Conduct the preparation according to the brightfield observation procedure.

4.3.2 Attach the dark condenser into the condenser holder the same way as Abbe condenser.(Ref, Page 8, 2.2.5)

4.3.3 Start your dark field observation.

4.4 Digital Camtra Head and Computer Observation(Fig.4)

4.4.1 Firmly screw the C-mount (SX2) into the digital camera head.

4.4.2 Loosen the clamping screw(B in Fig.3) on the trinocular observation tube H(A in Fig.3).

4.4.3 Fit the C-mount (SX2) into the Trinocular observation tube H and tighten the clamping screw.

4.4.4 Connect the digital camera head and the computer by USB cord.

4.5 Digital Camera Photography(Fig.4)

4.5.1 Firmly screw the Digital-Camera-Adapter(SD2+SDK) tightly.

4.5.2 Connect the digital camera with the SDK, then connect SD2 with SX2.

4.5.3 Loosen the clamping screw(B in Fig.3) on the trinocular observation tube H(A in Fig.3).

4.5.4 Fit the SX2 into the trinocular observation tube H and tighten the clamping screw.

4.6 Film Camera Photography(Fig.4)

4.6.1 Screw open the Photo-Adapter(SY2).

4.6.2 Insert the required photo eyepiece into SY2 and screw close the SY2 tightly.

4.6.3 Connect the camera with SY2.

4.6.4 Loosen the clamping screw(B in Fig.3) on the trinocular observation tube H(A in Fig.3).

4.6.5 Fit SY2 into the trinocular observation tube H and tighten the clamping screw.

5.TROUBLESHOOTING GUIDE

Under certain conditions, performance of the microscope may be adversely affected by factors other than defects. If problems occur, please review the following list and take remedial action as needed.

If you cannot solve the problem after checking the entire list, please contact your local UOP representative for assistance.

problem	causes	Remedy
1) The bulb does not light.	Power cord of the power supply unit is unplugged.	Plug in the power cord into a power outlet.
	Main switch of the power	Set the main switch to "I" (ON)

	supply unit is not ON.	
	The fuse is burnt out	Replace the fuse.
	The bulb is burnt out.	Replace the bulb.
2) The bulb lights but the field of view is dark.	The voltage is too low	Increase light intensity to an optimum voltage.
	Condenser is not well positioned.	Adjust the condenser height until the field iris diaphragm image is formed in the specimen plane.
	Condenser is not centered.	Center the condenser so that the field iris diaphragm image is centered in the field of view.
	Revolving nosepiece is not in a click position.	Make sure that the revolving nosepiece clicks properly into place.
	Field iris diaphragm is not opened wide enough.	Open the field iris diaphragm sufficiently.
	Too many filters are used.	Reduce the number of filters to the minimum required.
3) Field of view is obscured or not Evenly illuminated.	The objective that falls outside the condenser's illumination range is used.	Use a condenser that matches the objective.
	Field iris diaphragm is not properly centered.	Center the field iris diaphragm correctly.
	Field iris diaphragm is stopped down too far.	Open the field iris diaphragm sufficiently.
	Revolving nosepiece is in an intermediate position	Engage the revolving nosepiece at a click stop.
	A filter is stopped in an intermediate position.	Set the filter at the appropriate position.
	The frost filter is not engaged.	Engage the frost filter.
4) Dirt or dust is visible in the field of view.	Dirt/dust on the specimen.	Clean thoroughly.
	Dirt/dust on the eyepieces.	
	Dirt/dust on a mirror unit.	
	Dirt/dust on the optical element.	
	Condenser is not correctly positioned and the frosted filter or filter is focused.	Adjust the condenser height until the field iris diaphragm image is formed in the specimen plane.
5) Image glares	Condenser is raised too high.	Lower to the proper position.
	Aperture iris diaphragm is stopped down too far.	Open the aperture iris diaphragm.
6) Visibility of observe image is	Objective in use is not designed for UCIS series.	Replace with an objective designed for UCIS optics.

poor, Image is not sharp, Contrast is poor, Details are poorly visible.	Front lens of the objective is dirty	Clean the objective.
	Immersion oil is not being used with an oil immersion objective.	Use immersion oil.
	Inappropriate slide or cover glass thickness.	Replace with glass of appropriate thickness.
	Dirt/dust on glass components (condenser, objective, eyepieces,)	Clean thoroughly.
	Phase plate are not centered.	Center it.
7) Image is blurred.	Objective is engaged incorrectly in the light path.	Make sure that revolving nosepiece clicks into place correctly.
	Specimen is tilted with respect to the stage.	Place the specimen correctly on the stage.
8) Field of view of one eye does not match that of the other.	The interpupillary distance is incorrect.	Adjust the interpupillary distance.
	Incorrect diopter adjustment.	Adjust the diopter.
	You are not accustomed to parallel optical axis.	When looking into eyepieces, do not stare at image from the beginning but see the overall field of view. It is sometimes recommended to turn your eyes away from the eyepieces, look far off and look in to the eyepieces again.
9)The coarse/ine adjustment knobs will not rotate easily or at all.	The rotation tension adjustment ring is too tight.	Loosen the ring optimally.
10) The stage lowers by its own weight.	The rotation tension adjustment ring is too loose.	Tighten the ring optimally.

BAM200 Series Main Modules and Specifications

Optical System	UCIS Infinity Independent Achromatic Optical Design
Observation Tube	30° inclined, Interpupillary Distance 52–74mm, 360° rotatable, Trinocular light is split 50/50 by a high quality coated prism.
Eyepieces	10X, 21mm High eye point, View field ϕ 18mm, Diopter adjustable, Eye shade optional.
Nosepieces	Inward, Four Position.
Focusing system	Working distance 60mm, Low positioning, Co-axial coarse and fine focusing knob, Coarse knob tension adjustable, With focus stop, Coarse:14mm per rotation, Fine:0.1mm per rotation, Minimum reading:1 microns on left-side focusing Knob.

Mechanical Stage	(156mmx138mm) platform with X-Y travel of 76 mm x 54mm by low positioned X/Y coaxial control knob, with scale mark and specimen-slide clip.
Illumination	Built-in 6V/20W halogen bulb, Optional with field diaphragm provides, Koehler type illumination for clear, sharp images. The brightness intensity control is conveniently located and provides steady even illumination at all levels.
Electric	CE approved electrical parts, three wire grounded electric cord.
Condensers	Abbe condenser NA 1.25 with iris diaphragm and with socket to accomodate PH slide. Darkfield condensers is available and can be installed quickly and easily.

BAM200 Series Microscope Standard Outfits

Modules & Specification				BAM202	BAM203	Quantity
BAM200 Frame				●	●	1
Eyepieces	GM10X/20	Plan Eyepiece 10X FN20		●	●	2
	GM16X/14	Plan Eyepiece 16X FN14		○	○	2
	GPM10X/20	Eyepiece with Reticle		○	○	1
	WY1	Centering Telescope		○	○	1
Slide Binocular Tube	Inclined30, rotatable 360, interpupillary distance range 52–74mm.			●	○	1
Slide Trinocular Tube	Inclined30, rotatable 360, interpupillary distance range 52–74mm.			○	●	1
Objectives	Infinity E-Plan	4X		●	●	1
		10X		●	●	1
		40X	(Spring)	●	●	1
		100Xoil	(Spring)	●	●	1
	Infinity Plan	10X		○	○	1
		20X		○	○	1
		40X	(Spring)	○	○	1
		100Xoil	(Spring)	○	○	1
Phase Contrast Slide	Infinity E-Plan Phase Contrast	10X		○	○	1
		40X	(Spring)	○	○	1
		100Xoil	(Spring)	○	○	1
	Infinity Plan Phase Contrast	10X		○	○	1
		40X	(Spring)	○	○	1
		100Xoil	(Spring)	○	○	1
Condenser	AbbeN. A. 0.90/1.25oil. Center adjustable, with aperture iris diaphragm, height adjustable			●	●	1

Dark Condenser			○	○	1
Bulb	6V20W, Osram Halogen Bulb		●	●	1
Kohler with Field Diaphragm			●	●	1
C-Mount	SX2		○	○	1
Digital Camera Adapter	SD2+SDK		○	○	1
Photo Tube	SY2		○	○	1
Photo Eyepiece	S2.5	2.5X	○	○	1
	S4	4X	○	○	1
	S6.3	6.3X	○	○	1
Analyzer			○	○	1
Polarizer			○	○	1
Filter	Blue		●	●	1
	Yellow		○	○	1
	Green		○	○	1
Filter Holder (Attached to the condenser holder, movable)			●	●	1
Filter Holder (Attached to the Kohler)			○	○	1
Immerse Oil			●	●	1
Allen Screw Driver			●	●	1
Power Cord	UOP supplied meets local power requiremmt		●	●	1
Symbol meaning:● means Standard:○ means Optional					