

# USER MANUAL

## OPERATING LAMPS



|         |                                                  |
|---------|--------------------------------------------------|
| LD05.06 | OPERATING LAMP, MOBIL                            |
| LD10.01 | OPERATING ROOM CEILING LAMP - DIGITAL SINGLE     |
| LD15.11 | OPERATING ROOM CEILING LAMP - ANALOG SINGLE      |
| LD10.02 | OPERATING ROOM CEILING LAMP - DIGITAL DUAL       |
| LD15.22 | OPERATING ROOM CEILING LAMP - ANALOG DUAL        |
| LD10.03 | OPERATING ROOM CEILING LAMP - DIGITAL            |
| LD15.33 | OPERATING ROOM CEILING LAMP - ANALOG             |
| LD20.19 | OPERATING ROOM CEILING LAMP WITH 19" MONITOR     |
| LD20.24 | OPERATING ROOM CEILING LAMP WITH 24" MONITOR     |
| LD20.52 | OPERATING ROOM CEILING LAMP WITH MONITOR         |
| LD20.53 | OPERATING ROOM CEILING LAMP, ARM MOUNTED MONITOR |



**INSPIRAL**



## TABLE OF CONTENTS

|        |                                                   |       |
|--------|---------------------------------------------------|-------|
| 1.     | GENERAL INFORMATION .....                         | 4     |
| 2.     | INTRODUCTION .....                                | 5     |
| 3.     | SAFETY .....                                      | 6     |
| 3.1.   | Symbols And Explanations.....                     | 7     |
| 3.2.1. | Sample Label .....                                | 7     |
| 4.     | PACKAGING, HANDLING AND STORAGE.....              | 08-09 |
| 4.1.   | Physical Environment.....                         | 8     |
| 4.2.   | Handling Conditions.....                          | 9     |
| 5      | PRESENTATION.....                                 | 10-12 |
| 5.1.   | Intended Purpose And Places Of Usage .....        | 10    |
| 5.2.   | Device And Part Descriptions.....                 | 10-11 |
| 5.3.   | Technical Properties .....                        | 13    |
| 6.     | INSTALLATION AND USE.....                         | 14-29 |
| 6.1.   | Ceiling Assembly.....                             | 14    |
| 6.1.1. | Flange Assembly.....                              | 14    |
| 6.1.2. | Pipe Flange Assembly .....                        | 14    |
| 6.1.3. | Horizontal Rotary Arm Assembly .....              | 15    |
| 6.1.4. | Hood Assembly.....                                | 16    |
| 6.1.5. | Rod Connection System .....                       | 17    |
| 6.1.6. | Spring Arm Assembly.....                          | 18    |
| 6.1.7. | Semilunar (Stabilizer) Arm Assembly.....          | 19    |
| 6.1.8  | Monitor Assembly (Duo LED CC) .....               | 20    |
| 6.2.   | Control Panel Demonstration .....                 | 21-25 |
| 6.3.   | Braking Adjustment In Horizontal Plane.....       | 26    |
| 6.4.   | Braking Adjustment In Vertical Plane .....        | 27    |
| 6.5.   | Focus on Diameter Enlargement and Reduction ..... | 28-29 |
| 7.     | CLEANING.....                                     | 30    |
| 7.1.   | Reflector Glass And Control Panel Cleaning .....  | 30    |
| 7.2.   | Sterilization Of The Light Handle .....           | 30    |
| 8.     | MAINTENANCE.....                                  | 31    |
| 8.1.   | Led Ampoule Change .....                          | 32-33 |
| 8.2.   | Possible Failures And Remedies.....               | 34    |
| 8.3.   | Electric Circuit Diagram.....                     | 35    |
| 9.     | RECYCLING .....                                   | 36    |
| 10.    | AUTHORIZED SERVICES.....                          | 37    |

**Copyright**

This document was prepared by (INSPITAL Medical) and technical specifications of products may be modified. This document cannot be reproduced or copied in whole or part without prior written permission of the Company.

**CE Compliance Declaration**

Operating Lamps meet requirements set forth in “General Rules for Electrical Medical Devices” EN 60601-1, EN 60601-2-41, and Medical Devices Directive MDD/93/42/EEC.

**Notified Body Information**

Name :Kiwa Belgelendirme Hizmetleri A.Ş.

Adress :(İTOSB) İstanbul Tuzla Organize Sanayi

Bölgesi Tepeören Mevkii

34957 Tuzla-İstanbul Turkey

Phone :+90 216 593 25 75

Fax :+90 216 593 25 74

**Patent**

Trade name and mark belong to INSPITAL Inc.

**Warranty Period**

INSPITAL branded LED Lamps are under guarantee for 2 (two) years against material, workmanship and production faults including all the parts unless it is used as shown in the user guide. Below mentioned situations are out of the scope of guarantee:

1. Faults and damages caused by usage errors,
2. Faults and damages occurred during the loading, unloading and transportation after the delivery of the product to the customer,
3. Faults and damages that may arise due to the brown out or voltage rise; faulty electrical installation, using the device at a voltage level different from the written voltage values on the label of the product,
4. Faults arising from using the product non-compliantly to the matters specified in the user guide.

Thank you for purchasing INSPITAL branded Operating Lamp. Our product is rigorously designed and produced in order to provide you the best quality and performance. Please read all the information regarding operation and safety of this device cautiously and keep this manual as a future reference guide.

This guide describes the Installation, Operation, Cleaning and Maintenance steps for INSPITAL branded Operating Lamp.

## **Points to consider during transportation and handling**

Please consider the instructions written on the device packages during transportation and handling.

## **Matters regarding to maintenance and repair**

Please call the authorized device service for the cases requiring maintenance and repair.

## **Information regarding to connection and assembly**

Please consult to the consequent departments for connection and assembly information.

**In order to have the lights used by providing the safety of the patient and the user and to have their maintenance done:**

- Before taking the light into operation, be sure that this user guide is read and understood together with all the warnings and explanations.
- Buyer is responsible for ensuring that the personnel, who is responsible for use and maintenance of Pendant System, is informed of content of this user's manual.
- Keep this user guide at such a place that the personnel responsible for the operation and maintenance of the light can reach.
- Do not use the light except the purposes specified in this guide.
- Perform the maintenance of the light as shown in the figure in Part 8.1 Maintenance.
- INSPITAL R&D department reserves the right to change and/or delete the technical properties of the device and the information given in this user guide.



**CAUTION!**

- Do not operate the device before you read the user guide.
- The device should only be used by authorized personnel.
- Comply with the warnings and instructions given in the user guide. The device should only be used by authorized personnel.

### 3.1 SYMBOLS AND EXPLANATIONS

| Symbols                                                                             | Description                                                    | Symbols                                                                             | Description                        |
|-------------------------------------------------------------------------------------|----------------------------------------------------------------|-------------------------------------------------------------------------------------|------------------------------------|
|    | Type B device                                                  |    | Dangerous voltage                  |
|    | Caution! Consult to the accompanying documents.                |    | Protective grounding               |
|    | Casing can only be opened by trained service personnel         |    | Equipotentiality                   |
|    | AP Category Device according to IEC 601                        |    | Do not carry by holding from here. |
| IPX4                                                                                | Protection stages provided by the casings according to IEC 529 |    | Do not invert.                     |
|    | Maximum carriage capacity                                      |    | Do not invert.                     |
|  | Comply to MDD 93/42/EEC Medical Devices Ordinance              |  | Caution! Fragile                   |
|  | Date of production                                             |  | Ambient Temperature                |

### 3.2.1. SAMPLE LABEL



İNSPİTAL Medikal ve Teknoloji A.Ş.  
Karaoğlu Mah. Küme Evleri No:745 06830  
Gölbaşı / Ankara, TÜRKİYE  
T: 0.312.619 02 22 F: 0.312.619 02 25  
www.inspital.com.tr

**Ürün Adı: AMELİYATHANE TAVAN LAMBASI ANALOG - İKİLİ**  
**Product Name: OPERATING ROOM CEILING LAMP - ANALOG DOUBLE**  
 Model : LD15.22  
 Giriş Voltajı (Input Voltage): 220-240 VAC / 50-60 Hz Güç (Power): 2x75W  
 Çalışma Voltajı (Operating Voltage): 24 VDC

**SN** 00000000000   
 8682079504324

2018      



#### CAUTION!

No amendments should be done on the device label and the label should not be re moved.

### 4.1. PHYSICAL ENVIRONMENT

Ambient temperature: 0°C ~ +50°C

Humidity: %40

**Operating lamp is dispatched by supporting with soft materials inside the package. After the light is delivered;**

- Carefully open the closed box from its appropriate locations.
- Carefully remove the nylon packaging.
- Carefully inspect the operating lamp and its accessories and determine whether they are damaged during transportation or not. If they are damaged during transportation, consult to the relevant sales point for your operating lamp and/or accessories.
- Carefully remove the accessories from their packaging.
- Store the operating lamp inside a box only by adding up maximum two boxes.



#### CAUTION!

The direct contact of other devices (aspirator, cautery ,etc) that are present in the environment where the operating room light is used with the light casing. Otherwise it may damage the light casing.

#### 4.2.HANDLING CONDITIONS

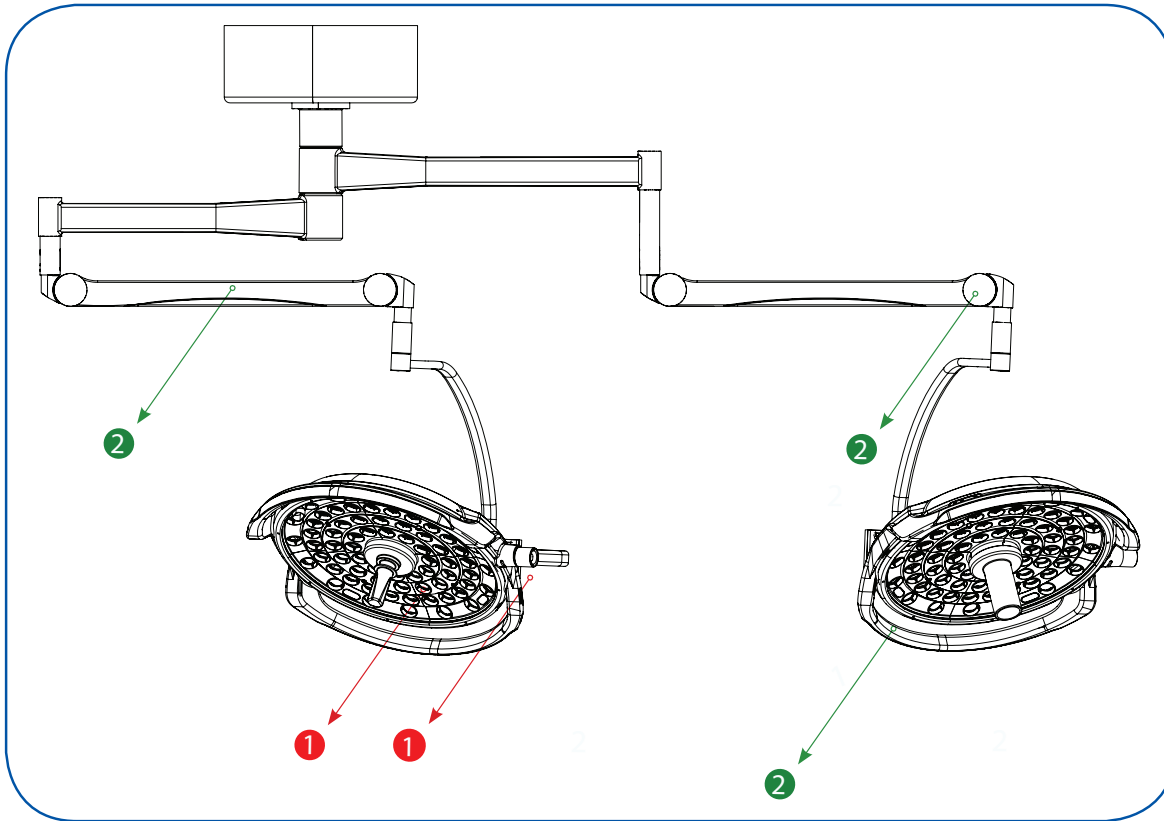


Figure 1



#### CAUTION!

While opening the nylon packaging of the operating lamp; never use sharp objects such as knives etc. otherwise you may damage the operating lamp.


Never carry the operating lamp by supporting from no. 1 locations shown in figure 1. In order to carry the operating lamp by supporting from no. 2 locations shown in figure 1.



## 5.1. INTENDED PURPOSE and PLACES OF USAGE

### INSPIRAL branded operating lamps;


- provide to view a real-like a detailed operation area,
- provide bright, shadeless, reliable and high-quality lightening due to its versatile special reflector system,
- provide flexible movement capability and perfect positioning opportunities due to its arms connected to a single center and its rack,
- provide to monitor the educational surgeries and to record and document them by its camera and monitor options,
- The light armature is designed to provide opportunity for the movement of the light cover and arranged in order to stay fixed where the light cover is left.

## 5.2. DEVICE and PART DESCRIPTIONS

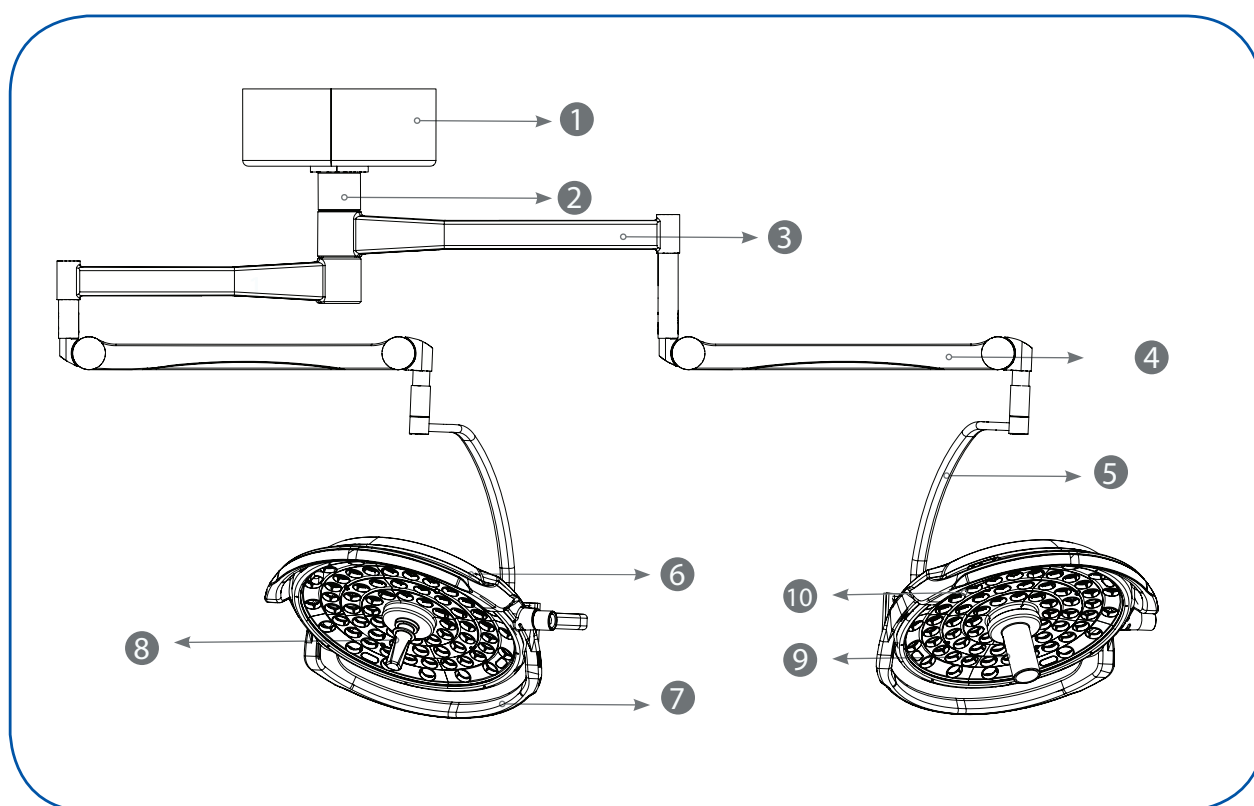

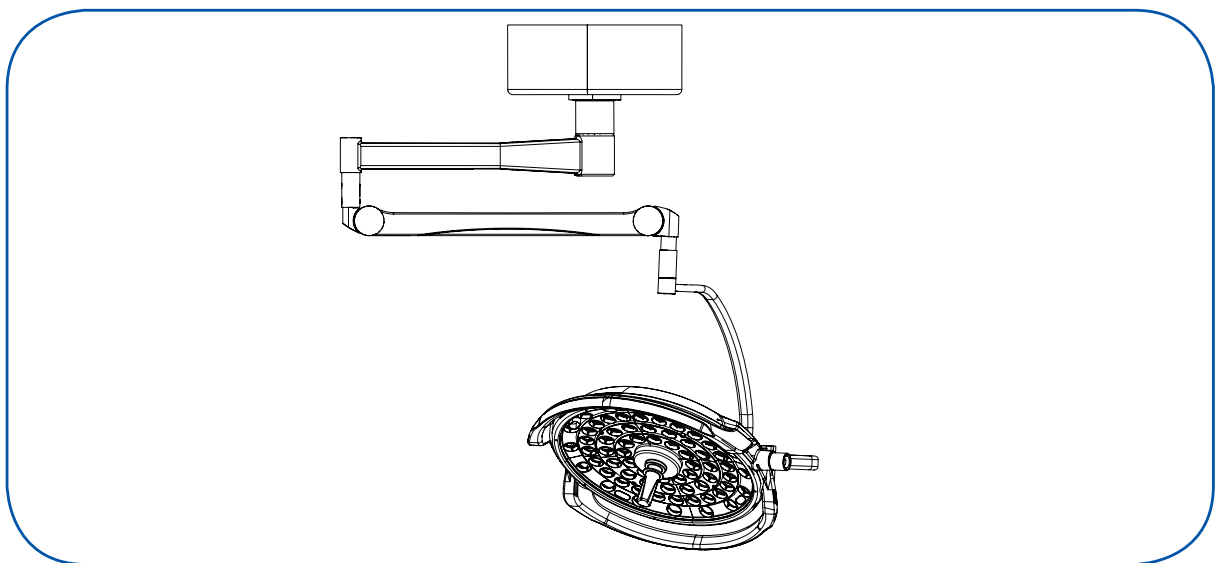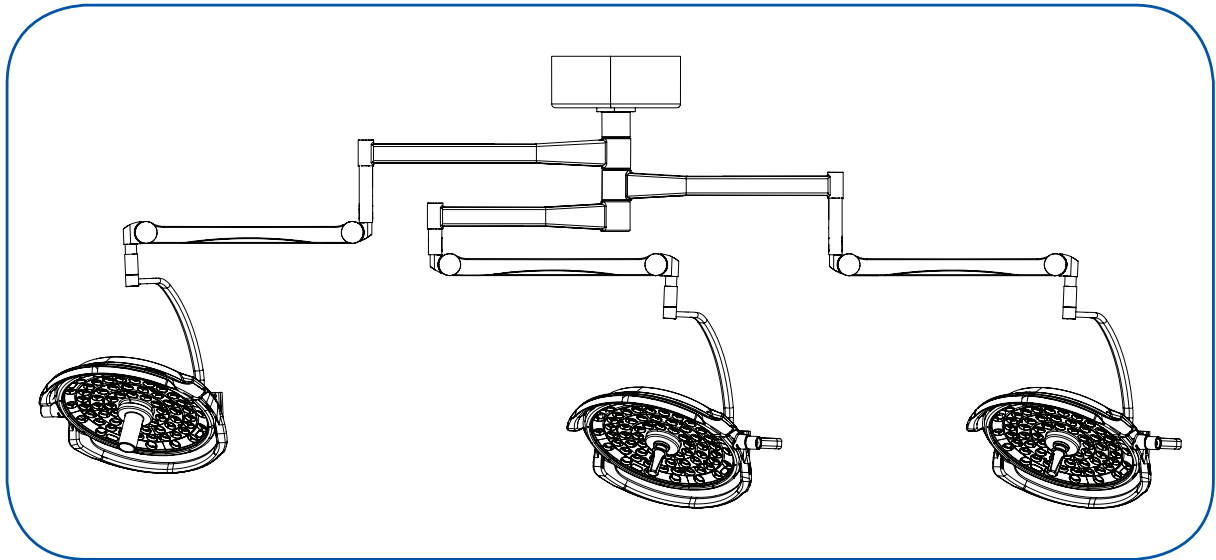

| OPERATING LAMPS |                            |         |                   |
|-----------------|----------------------------|---------|-------------------|
| PART NO         | PART NAME                  | PART NO | PART NAME         |
| 1               | HOOD                       | 6       | CONTROL PANEL     |
| 2               | FIXED CONVEYOR             | 7       | LIGHT COVER       |
| 3               | HORIZONTAL ROTARY ARM      | 8       | LIGHT HANDLE      |
| 4               | SPRING LEVER               | 9       | CAMERA (OPTIONAL) |
| 5               | SEMILUNAR (STABILIZER) ARM | 10      | LEDS              |

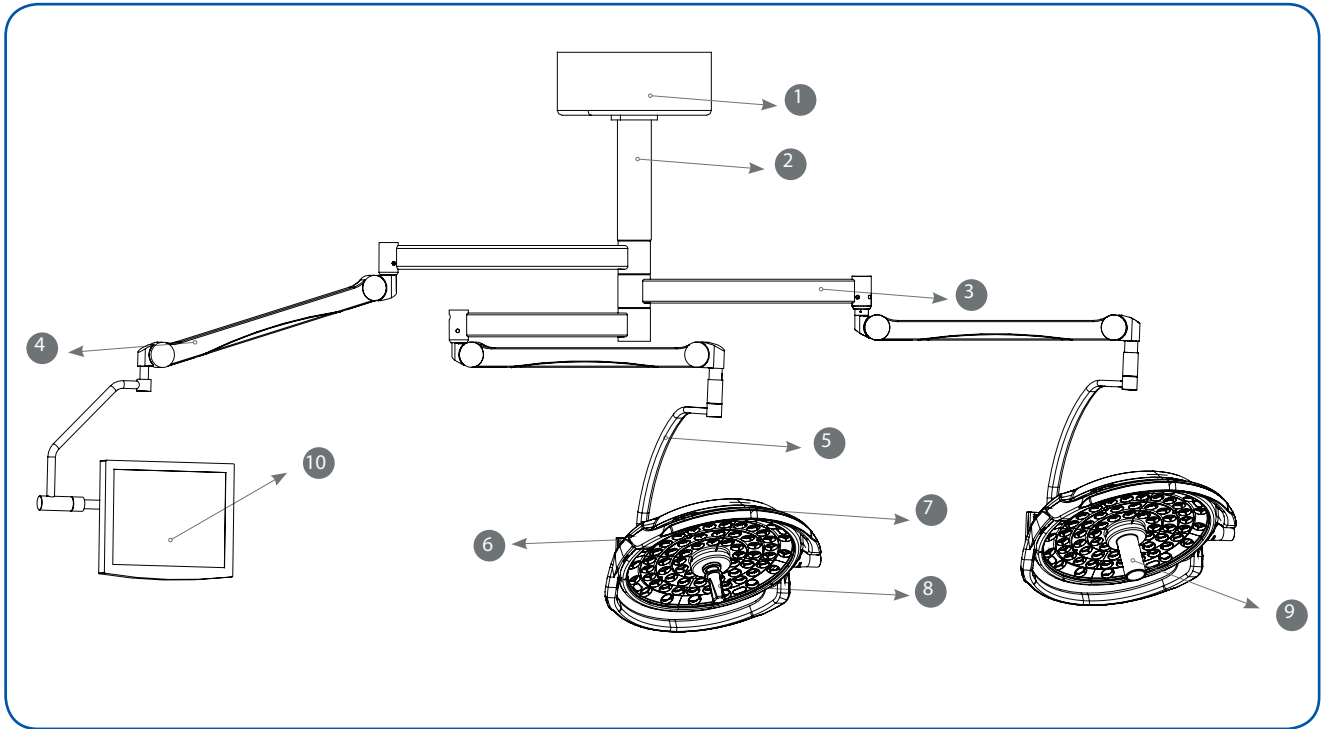

| OPERATING LAMPS |                         |         |                        |
|-----------------|-------------------------|---------|------------------------|
| PART NO         | PART NAME               | PART NO | PART NAME              |
| 1               | HOOD                    | 6       | CONTROL PANEL          |
| 2               | FIXED SUPPORT           | 7       | LAMP HEADGEAR          |
| 3               | HORIZONTAL ROTARY ARM   | 8       | LAMP (FOCUSING) HANDLE |
| 4               | SPRING ARM              | 9       | CAMERA (OPTIONAL)      |
| 5               | HALF MOON (BALANCE) ARM | 10      | MONITOR (OPTIONAL)     |



**DİKKAT!**

No. 9 and 10 camera and monitor connections are performed optionally.


### 5.3 TECHNICAL PROPERTIES

| Titles                                    | Main Lighthouse |
|-------------------------------------------|-----------------|
| Lightening cap                            | 600 mm          |
| Light power, 1 m                          | 160.000 lux     |
| Color separation index, CRI               | >96             |
| Light intensity adjustment                | 1 – 100 %       |
| Temperature increase in the surgical area | < 0,5 0°C       |
| LED module quantity                       | 68              |
| LED life                                  | >50.000 hours   |
| Total consumption                         | 85 W            |
| Lightening depth (light beam)             | 1300 mm         |
| UV-radiation level                        | 0,092 W/m2      |

| Monitor Specifications | 19" Monitör           | 24" Monitör           |
|------------------------|-----------------------|-----------------------|
| LCD Panels :           | Aktif Matriks TFT LCD | Aktif Matriks TFT LCD |
| Viewing area           | 19"                   | 24"                   |
| Resolution             | 1280 x 1024           | 1920 x 1200           |
| Display ratio          | 05:04                 | 16:10                 |
| Viewing angle          | 1700                  | 178                   |
| Response               | 14 ms                 | 14 ms                 |

| Camera Specifications            |
|----------------------------------|
| SONY HD Analog Camera            |
| 1/2.8 type CMOS sensor           |
| Full HD z1080p                   |
| 30x optik zoom                   |
| Digital / Analog HD video output |
| Auto Focus                       |

## 6.1.6 SPRING ARM ASSEMBLY

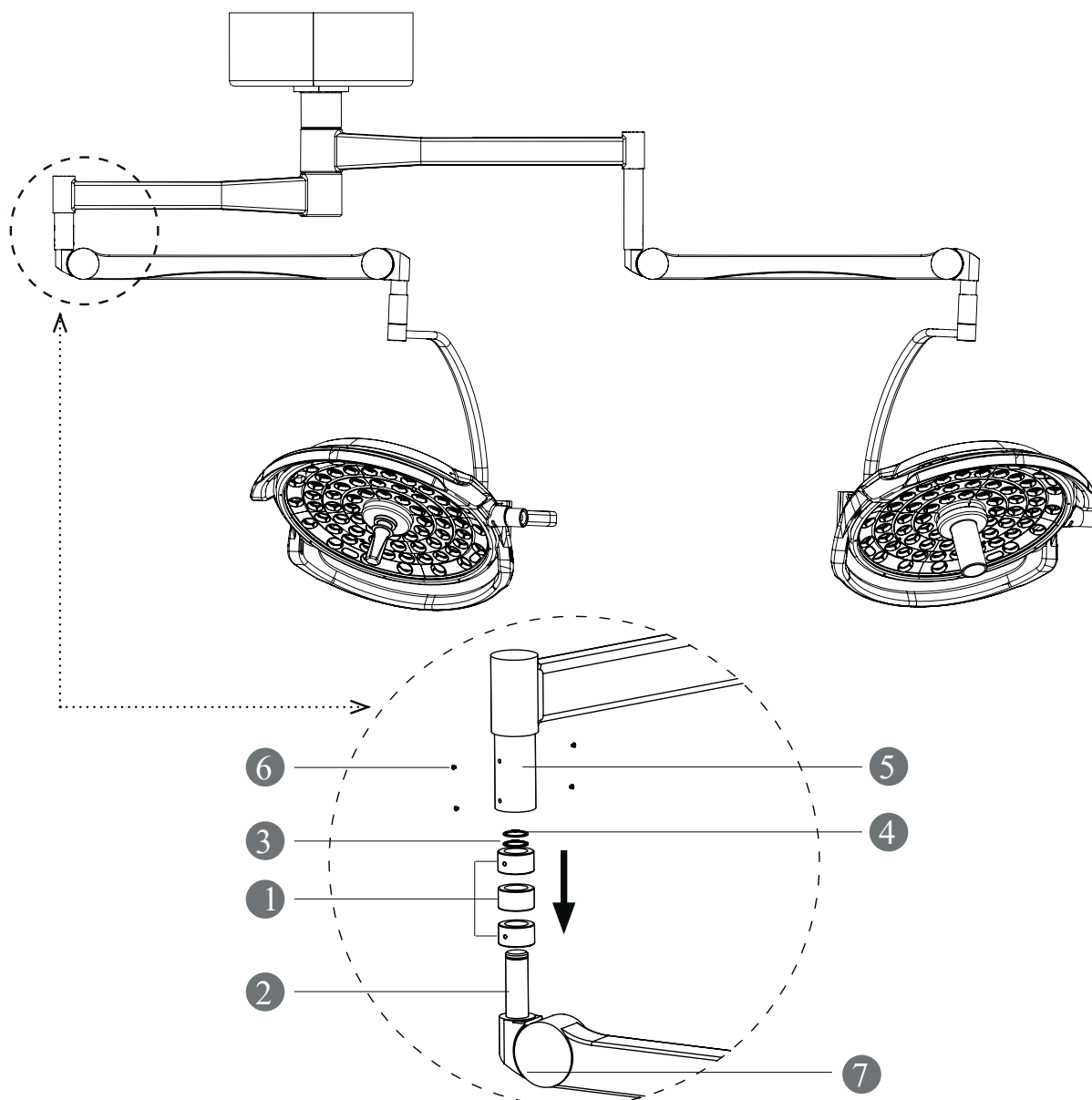

Fix the 3 set of clamps (1) from the light package to the shaft (2) at the end of the spring arm. Affix the washer (3) and the safety ring (4). Place the horizontal rotary arm (5) in the clamps and affix them with 4 screws (6).


At the end of this process; the electrical connections will be done with the sockets inside the shafts. In order to adjust the minimum and maximum movements of the ceiling height and the light, use No. 7 braking pin.

**CAUTION!**

Loosening, removing or misusing the braking pin (7) by the user may cause the complete elimination of the braking mechanism of the movable conveyor arms. The users should not interpose with the conveyor arms. In case of a fault; technical service should be called.

## 6.1.7 6. INSTALLATION AND USE

## 6. INSTALLATION AND USE

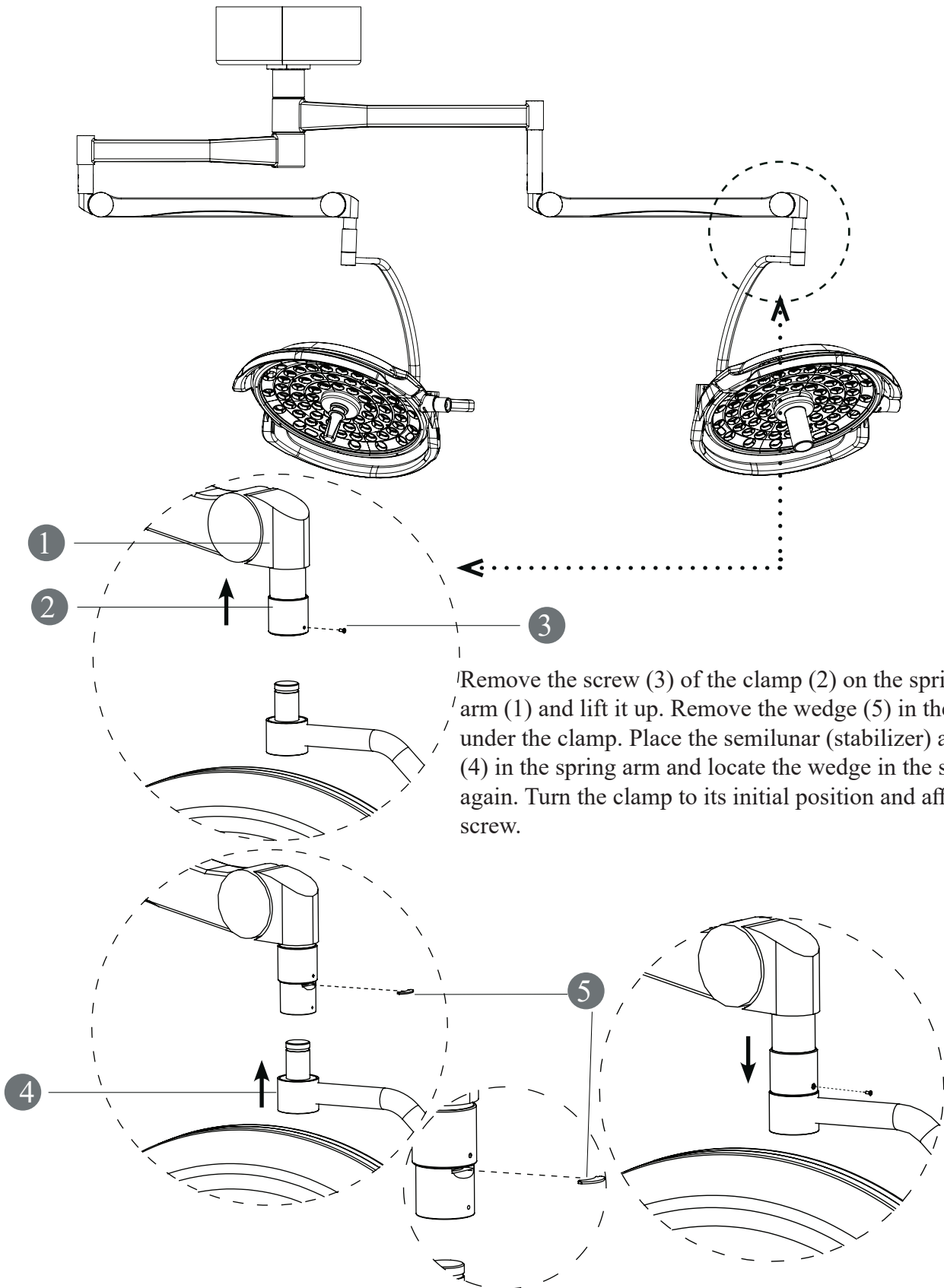


Remove the screw (3) of the clamp (2) on the spring arm (1) and lift it up. Remove the wedge (5) in the slot under the clamp. Place the semilunar (stabilizer) arm (4) in the spring arm and locate the wedge in the slot again. Turn the clamp to its initial position and affix its screw.

### WARNING!

Removal of no. 5 wedge from its location may cause vibrational impacts or have the existing suspended mass to fall. In case the wedge is loosened or removed; notify the technical service. Provide the safety of the patient and the user.

## 6.1.8 (DUO LED CC) MONITOR ASSEMBLY

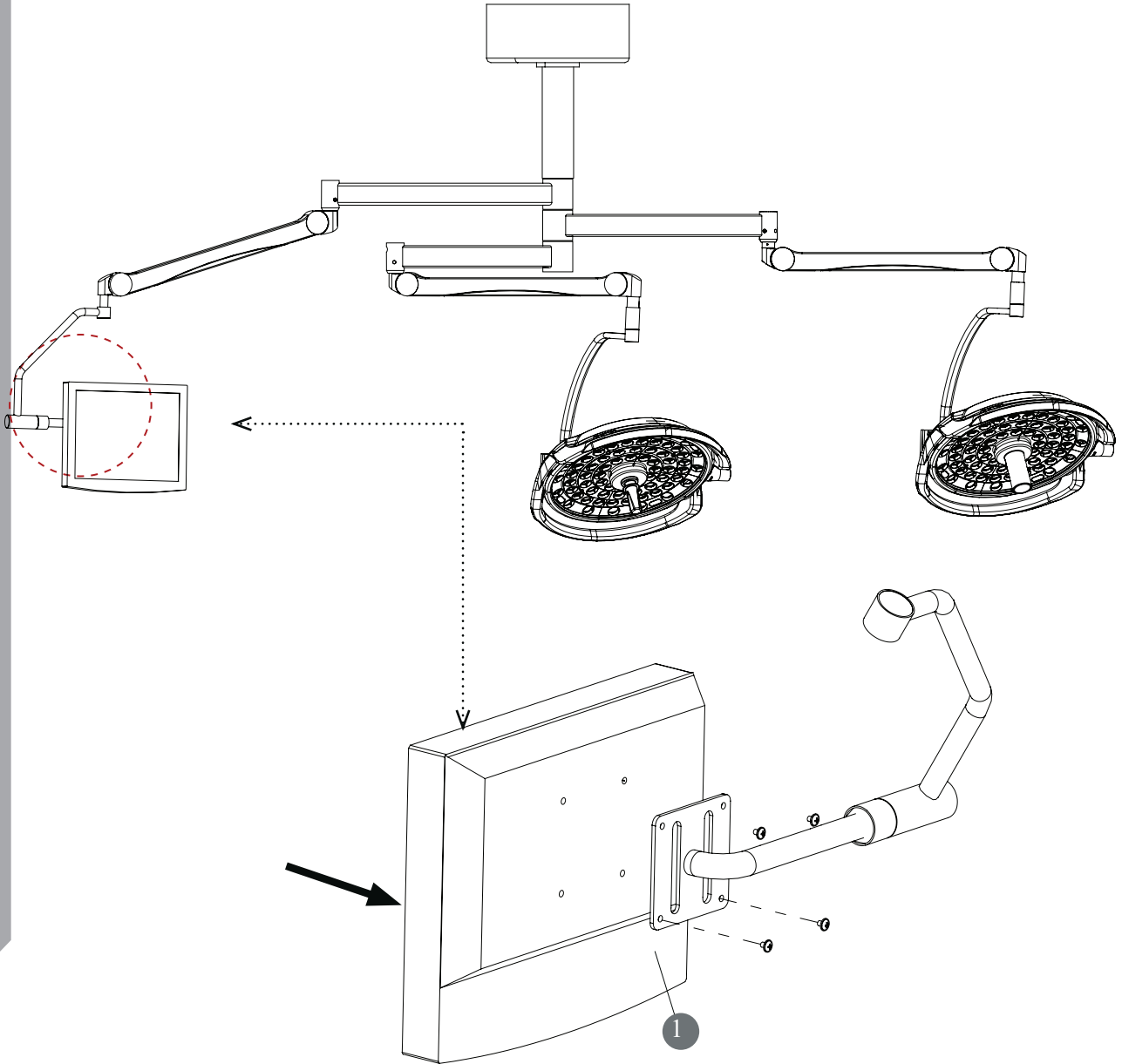

Mount the monitor to the arm flange and connect it with screws. Connect the cables that come out of the arm to appropriate monitor outlets. Connect the upper arm of the monitor as the arms of other lamp headgear.



**DİKKAT**


Monitor no. 1 is supplied optionally and should be assembled as shown.

6.2. CONTROL PANEL DESCRIPTION

SYSTEM OFF



SPLASH SCREEN

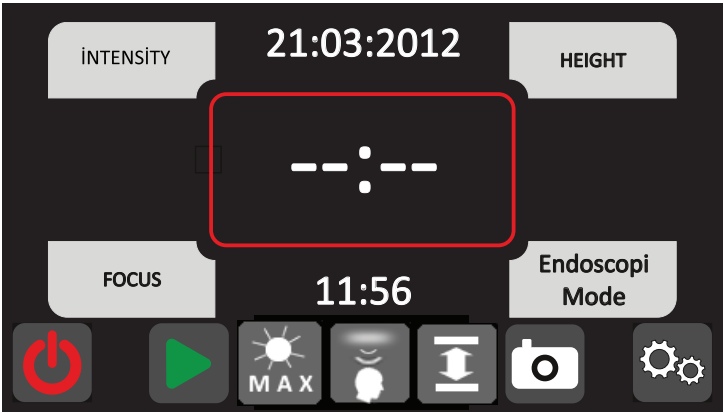

CLICKED ON START KEY



# 6. INSTALLATION AND USE


LIGHT POWER SETTING



FOCUS SETTING



ENDOSCOPY MODE

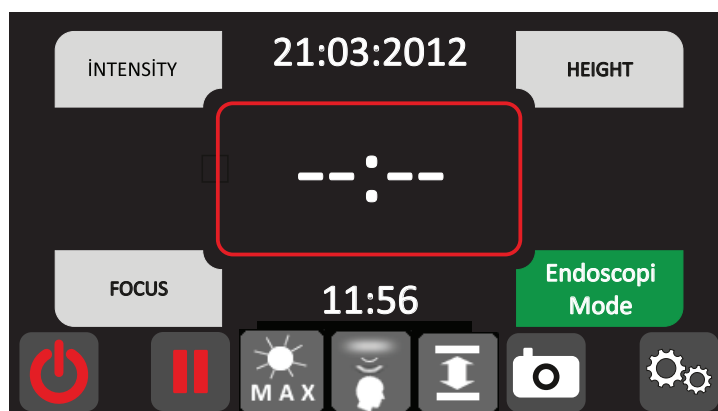

HEAD SENSOR OFF




↑  
FOR AUTOMATIC LIGHT SETTING  
ACC. TO THE HEIGHT


CAMERA ON




SERVICE SCREEN

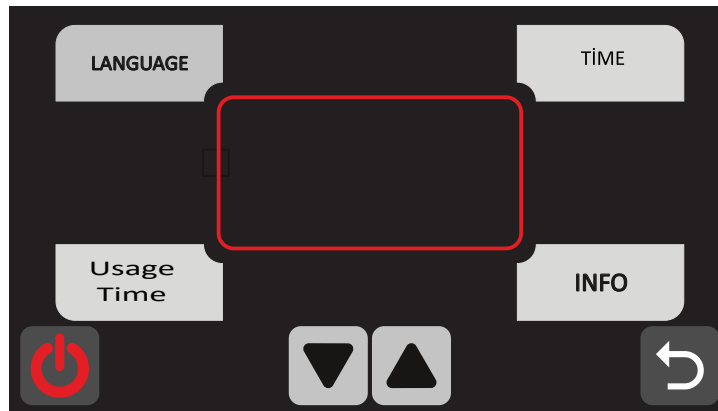



LANGUAGE OPTION

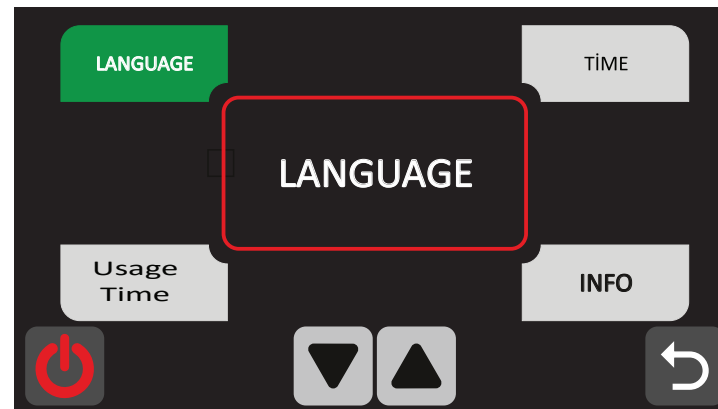


TIME SETTING

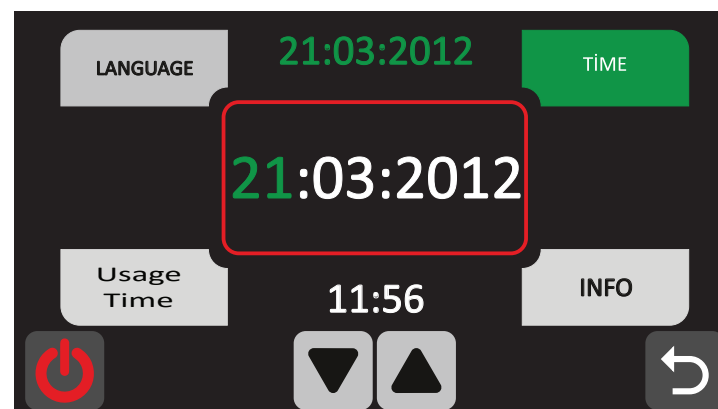

LIFETIME RESETTING

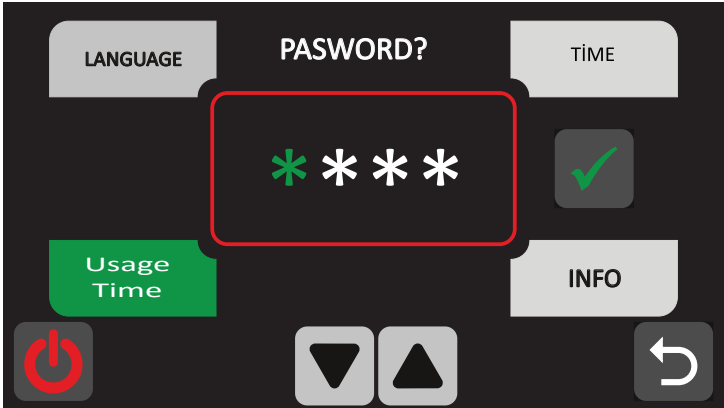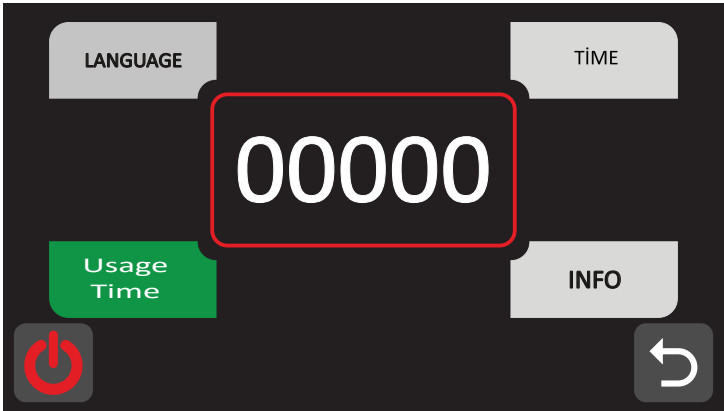

INFO SCREEN

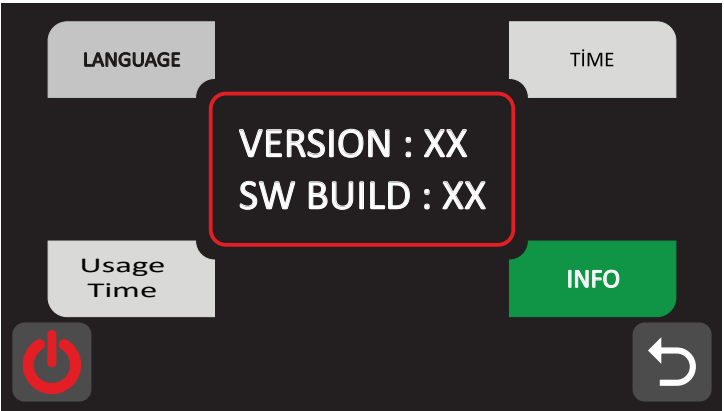


HEAD SENSOR ON




↑  
FOR SINGLE KEY 100% LIGHT POWER

## 6.3. INTRODUCTION OF ANALOGUE CONTROL PANEL

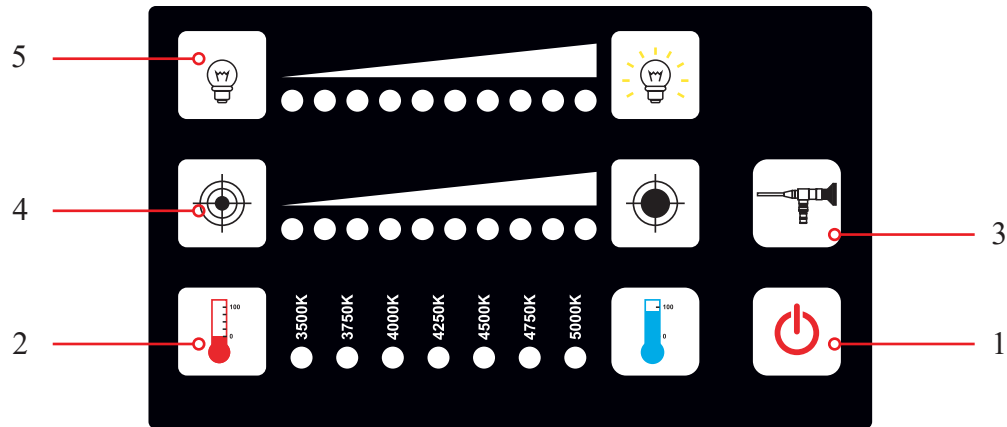


| NO | PARAMETER          | FUNCTION                                                            |
|----|--------------------|---------------------------------------------------------------------|
| 1  | On-Off             | Used for turning on and off the lamp                                |
| 2  | Color temperature  | To adjust color temperature between 3500-5000K arasında değiştirir. |
| 3  | Endoscopy          | To enable endoscopy mode                                            |
| 4  | Focus              | Used for adjusting focus diameter                                   |
| 5  | Illumination level | To adjust illumination level                                        |

## 6.3 BRAKING ADJUSTMENT IN HORIZONTAL PLANE

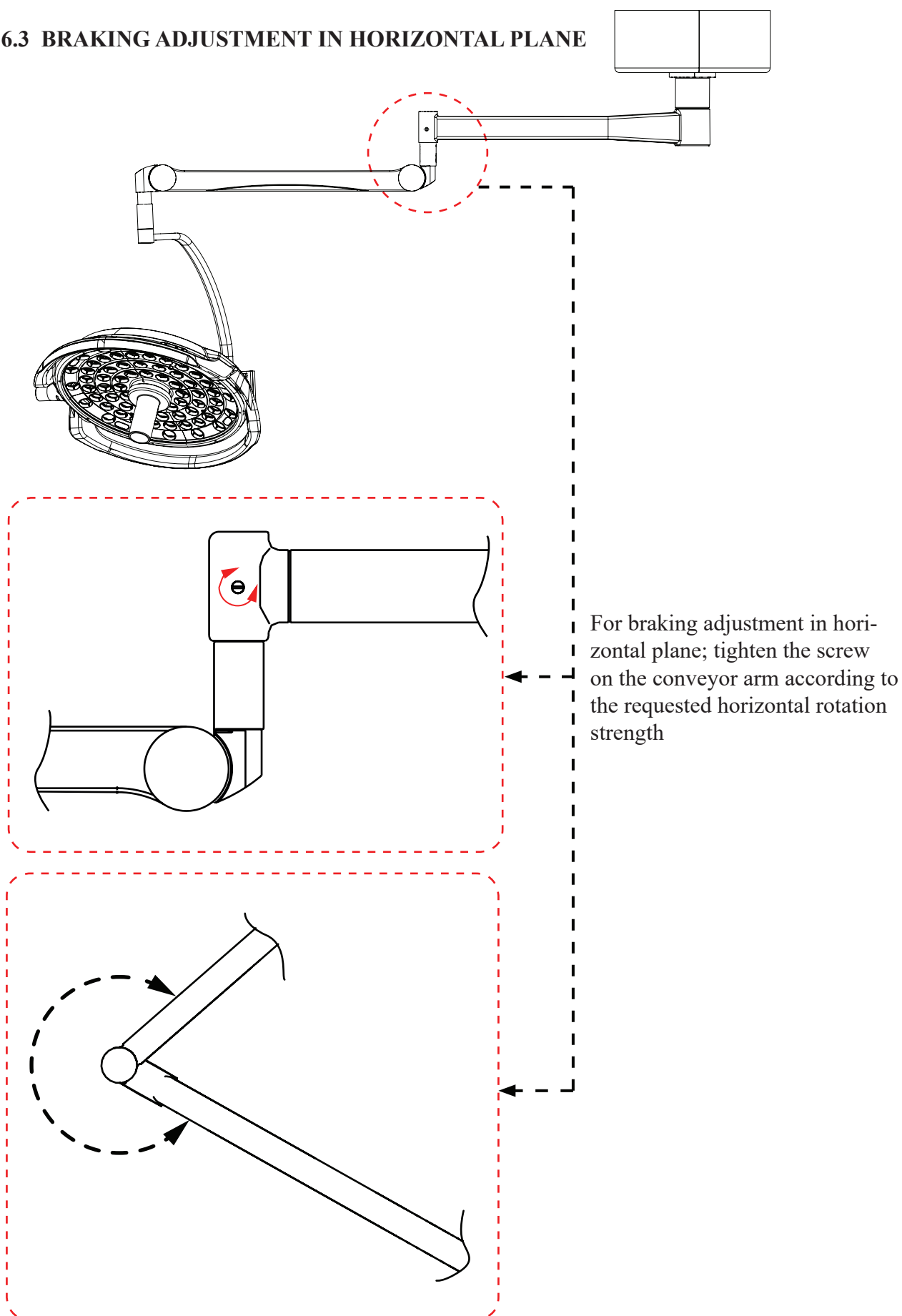

## 6.4 BRAKING ADJUSTMENT IN VERTICAL PLANE

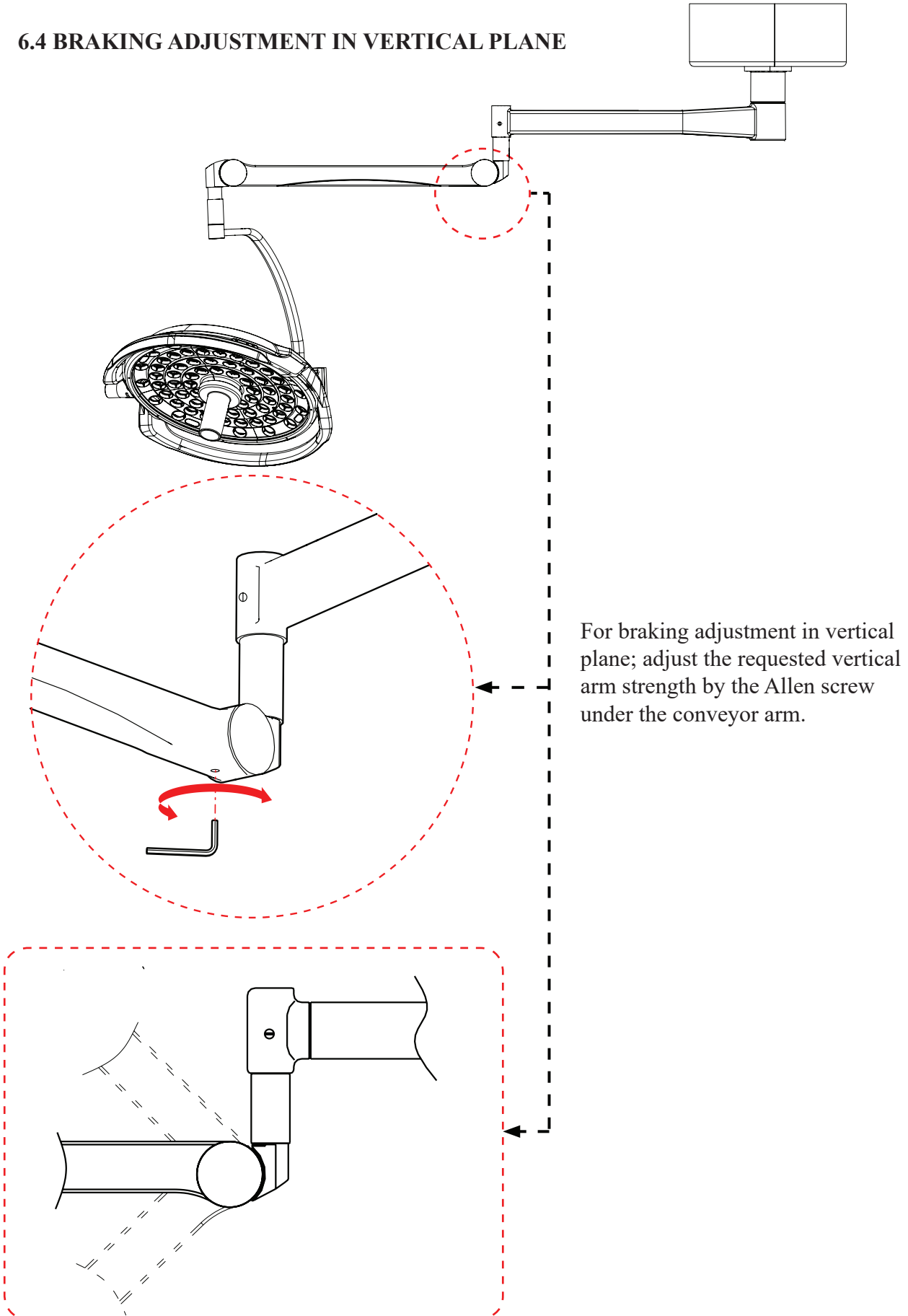

**CAUTION!**


For the operating room lights; the led ampoules with the properties specified in “Technical Properties” section of this guide. Using ampoules with different properties may cause temperature rises and may harm the patient and the user. In case there is any cracks, fractures and holes on the transparent plastic of the light lightening cap; the internal part of the light will become eligible for microorganism reproduction. Thus in case of the cracks, fractures and holes on the external side of the light; the material should be replaced with a new one.

**8.3. POSSIBLE FAILURES AND REMEDIES**

Please contact to our authorized technical service for other information that may help you for the operating lamp components and their repairs.

| Failure                                                      | Causes                                         | Remedies                                                                                                                                                                                                                                                                                                                                                     |
|--------------------------------------------------------------|------------------------------------------------|--------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|
| 1- Indicator on the control panel does not lit               | 1- Indicator on the control panel does not lit | <ul style="list-style-type: none"> <li>•Turn on the switch on the operation panel</li> <li>•Control the transformer warning light whether it is on or off. If off, replace the transformer.</li> <li>•If the warning light is on, control the cables on the arms.</li> <li>•In case the failure is not solved, consult to the authorized service.</li> </ul> |
| 2- Lamp couldn't be fixed and changes its position by itself | 2- The brakes might be loose.                  | The brakes must be tightened and control.                                                                                                                                                                                                                                                                                                                    |
| 3- Control box or operation panel defective                  | 3- Control box or operation panel defective    | Inform our company or Regional Failure Service.                                                                                                                                                                                                                                                                                                              |

8.4 ELECTRIC CIRCUIT DIAGRAM

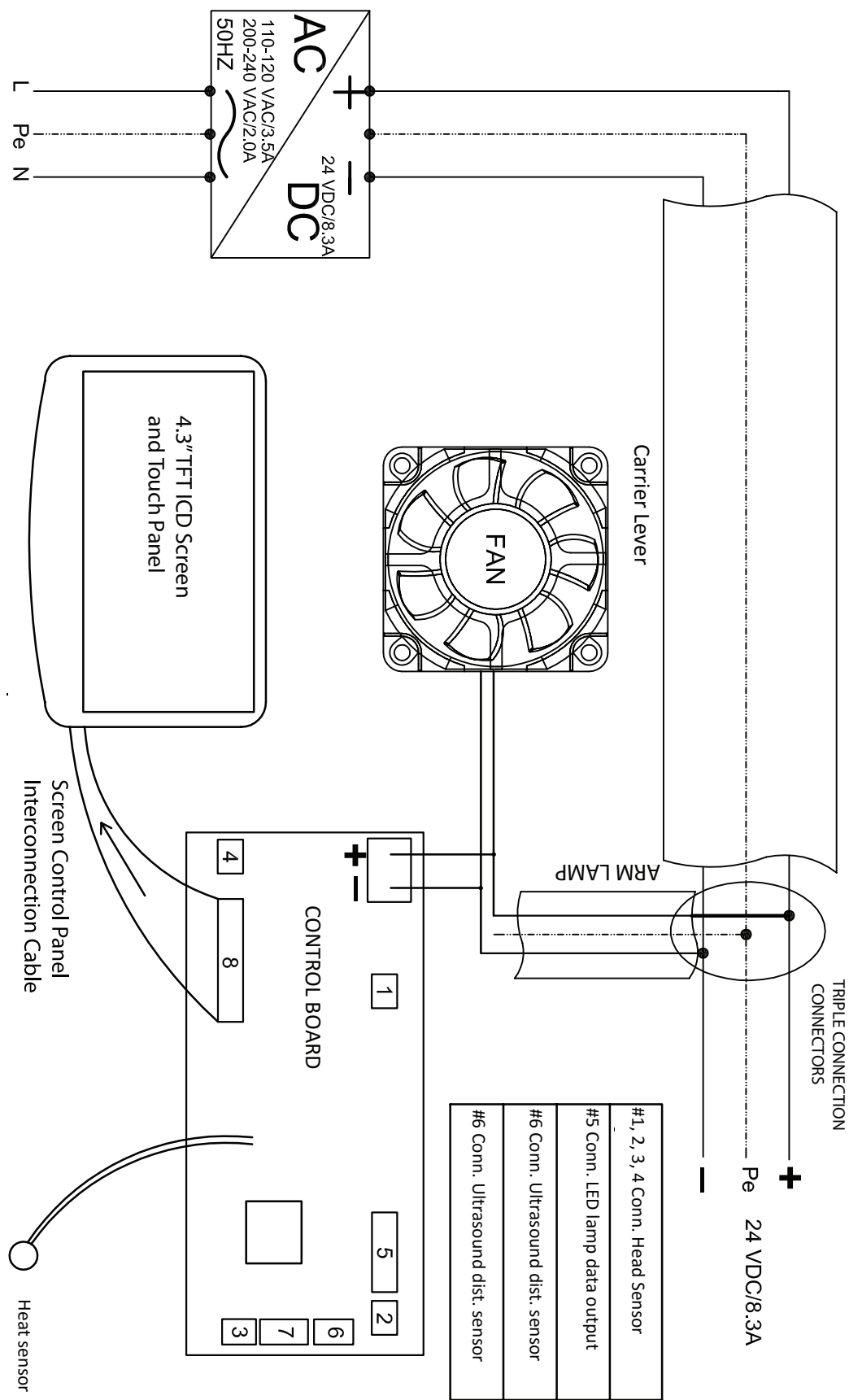



## 9. RECYCLING

If you will dispose or replace your operating lamp, then check recyclability of each part of it. Accessory frame is predominantly made of ABS plastic and aluminum. Both are recyclable materials.

In order to obtain further information about recycling please consult to relevant companies and plants or visit the sites on internet that give information on recycling.

## 10. AUTHORIZED SERVICES

**Ankara – Headquarters / Production Plant**

İnsipital Medikal Teknolojileri A.Ş

Karaoğlan Mahallesi Küme Evleri

No: 745 Gölbaşı / ANKARA/TURKİYE

Tel: +90 312 619 02 22

Fax: +90 312 619 02 25

[www.inspital.com..tr](http://www.inspital.com..tr)

# INSPI TAL

MEDICAL GAS SYSTEMS

---

## Headquarters & Factory

Karaoğlan Mahallesi Küme Evleri  
No: 745 Gölbaşı / ANKARA, TÜRKİYE  
Tel: +90 312 619 02 22  
Fax: +90 312 619 02 25

---

## İstanbul Region Headquarters

Oruç Reis Mah. Tekstil Kent Cad.  
Tekstil Kent Ticaret Merkezi  
10-AN NO: Z09 (A7 BLOK NO: 19)  
Esenler-İSTANBUL, TÜRKİYE  
[www.inspital.com.tr](http://www.inspital.com.tr)  
[info@inspital.com](mailto:info@inspital.com)