

Digitally signed by Belan Irina
Date: 2025.12.04 11:11:39 EET
Reason: MoldSign Signature
Location: Moldova

MOLDOVA EUROPEANĂ



SIEMENS
Healthineers

ACUSON Juniper
ACUSON Juniper Select
Diagnostic Ultrasound System
Instructions for Use

US



An updated User and Reference Manual may be available for your ultrasound system. Access the most recent manual from the Siemens Healthineers Document Library at doclib.siemens-healthineers.com.

Product Version 3.0

Software Version VC10

©2018-2025 Siemens Medical Solutions USA, Inc.

All Rights Reserved.

Date of first issue: 2025-08

All trademarks are the property of their respective owners.

Siemens Medical Solutions USA, Inc. reserves the right to change its products and services at any time. In addition, this publication is subject to change without notice.

Manufacturer

Siemens Medical Solutions USA, Inc.
Ultrasound
22010 S.E. 51st Street
Issaquah, WA 98029
U.S.A.
Phone: +1-888-826-9702
siemens-healthineers.com

Siemens Healthineers Headquarters

Siemens Healthineers AG
Siemensstr. 3
91301 Forchheim
Germany
Phone: +49 9191 18-0
siemens-healthineers.com

Contents

Chapter 1	Introduction Overview of the diagnostic ultrasound imaging system, intended use, and a detailed description of the touch screen and control panel.
Chapter 2	Safety and Care Detailed information on system safety and care of the ultrasound system. Includes information on system maintenance.
Chapter 3	System Setup Detailed descriptions for setting up and transporting the system. Includes procedures for updating your system software.
Chapter 4	Examination Fundamentals Procedures for registering a patient; selecting a transducer, exam type and operating mode; and optimizing and capturing images. Includes procedures for creating and editing user-defined exam types.
Chapter 5	Examination Review Procedures for reviewing images and clips in the current or previous examination.
Chapter 6	Transducer Care and Accessories Detailed information on transducer safety, cleaning and disinfecting, and transducer accessories. Attachment procedures for transducer accessories.
Chapter 7	Specialty Transducers Descriptions of the following specialty transducers: ▪ 9VE3 ▪ 10BCC3

Appendix A Technical Description

Technical description of the ultrasound system, including the operating and storage temperatures for the transducers and system.

Appendix B Patient Data

General description of the patient registration form and the patient browser.

Appendix C Cybersecurity

Information about ultrasound system security.

Appendix D Acoustic Output Reference

Acoustic output reporting tables.

Note: Not all features and options described in this publication are available to all users. Please check with your Siemens Healthineers representative to determine the current availability of features and options.

See also: For a list of available options and features, refer to Appendix A in the Instructions for Use.

About the User and Reference Manuals

The user and reference manuals contain descriptions for the following ultrasound systems:

- ACUSON Juniper diagnostic ultrasound system
- ACUSON Juniper Select diagnostic ultrasound system

Features and options unique to an ultrasound system are identified within the following locations in the user and reference manuals:

- Chapter 1, Instructions for Use
- Chapter 3, Instructions for Use
- Appendix A, Instructions for Use

The user and reference manuals consist of the following publications.

Publication	Includes
<i>Instructions for Use</i>	▪ Conventions and typographical conventions used in the manuals ▪ Intended audience ▪ Introduction to the ultrasound system ▪ Technical description of the ultrasound system ▪ Safety and care information for the system and compatible transducers ▪ Descriptions of system controls ▪ Procedures for system setup, examination fundamentals, image acquisition and optimization, and examination review ▪ Information about system security and protection of patient data against security threats ▪ Acoustic output data
<i>Clinical Applications Guide</i>	▪ Explanation of clinical applications and clinical software programs for use on the ultrasound system ▪ Procedures for the physiologic function and the biopsy function ▪ Description of the transesophageal transducer, including cleaning and care of the transducer ▪ General and exam-specific measurements and calculations ▪ Exam-specific annotations
<i>System Reference*</i>	▪ General and exam-specific measurement and calculation labels ▪ Description of customizable system settings ▪ Information about DICOM connectivity and network capabilities ▪ Clinical references ▪ Information regarding the electromagnetic compatibility (EMC) testing of this system

*Languages supported by the user interface include a translation of this publication.

Conventions

Take a moment to familiarize yourself with the following conventions in your user manual.

Convention	Description
Warnings and Cautions	 WARNING: Warnings are intended to alert you to the importance of following the correct operating procedures where risk of injury to the patient or system user exists.  Caution: Cautions are intended to alert you to the importance of following the correct operating procedures to prevent the risk of damage to the ultrasound system.
Notes and Prerequisites	Prerequisite: Prerequisites contain tasks the user must complete or information the user needs prior to performing a procedure. Note: Notes contain information concerning the proper use of the ultrasound system or correct execution of a procedure.
Examples of Cross-References	Examples: See also: For information about configuring a network, refer to Chapter 2 in the System Reference. See also: Touch screen controls, page 1-17 See also: For a list of transducers compatible with the needle guide accessories, refer to Chapter 6 in this manual.

Typographical Conventions

The following typographical conventions are in descriptions and procedures within your user and reference manuals. Use the conventions to identify the location of a control or selection on the ultrasound system.

Control or Selection	Description
Controls on the Control Panel and Keys on the Keyboard	Controls on the control panel are indicated by uppercase, boldface type. The function of the trackball and trackball keys is indicated by on-screen labels displayed on the image screen. Examples: ▪ Push DEPTH / FOCUS. ▪ Rotate 2D to adjust the 2D gain. ▪ Press CALIPER to activate the measurement function. ▪ Press and hold the Alt key on the keyboard to access the special characters. ▪ Press UPDATE. ▪ Position the marker and then press SET.
Controls on the Control Panel and Soft Key Selections on the Touch Screen	The touch screen has unlabeled controls located on the control panel. Tap, press, rotate or press and then rotate the control to make a selection or change a setting. The selection assigned to these controls is indicated by boldface type and enclosed in square brackets ([]). Examples: ▪ Press [Tx Power] downward to decrease transmit power. ▪ Rotate [Map].
Selections on the Image Screen	On-screen selections are indicated by boldface type. Roll the trackball to position the pointer on the on-screen control or object and then press a trackball key. Examples: ▪ Click Cancel. ▪ Double-click a thumbnail. ▪ Select the required check box to search for scheduled procedures.
Selections on the Touch Screen	On-screen selections are indicated by boldface type. To access the selection, tap the touch screen. To change the direction of an object displayed on the touch screen, drag the object on the touch screen. Examples: ▪ Tap Stress Echo. ▪ Drag the blue dot to indicate the direction of the arrow. ▪ Swipe the page to select the next page.


1 Introduction

System Overview	2
Intended Use	7
Indications for Use Statement.....	8
Operating Modes.....	9
Patient Target Groups	9
Intended Use Environment.....	9
Intended Users	10
Contraindications.....	10
Transducers and Intended Applications	11
Curved and Linear Array Transducers	11
Phased Array Transducers.....	16
Endocavity Transducers	17
<i>fourSight</i> 4D Transducers.....	18
CW Transducers	19
Legend	19
Principles of Operation	20
Clinical Benefits.....	20
Image Screen Layout.....	21
Screen Saver.....	21
System Status Indicators.....	22
About the Control Panel.....	24
Control Panel Lighting	25
Power Control.....	25
Workflow Controls	26
Operating Mode Controls.....	27
Image Acquisition and Documentation Controls	28
Optimization Controls	29
Retractable Keyboard.....	30
Special Keys.....	30
Special Characters	31
Touch Screen Controls.....	32
Touch Screen Symbols	33
Touch Screen Keyboard.....	34
Special Keys and Special Characters	34
System Configuration	35

System Overview


The ACUSON Juniper and ACUSON Juniper Select ultrasound systems are maneuverable, digital diagnostic ultrasound imaging systems. The systems utilize advanced imaging processing and transducer technology. The operating system is based on Windows technology. The system software supports standard applications, exam-specific imaging presets, measurements, body markers, annotations, reports, and system diagnostics.

See also: For a technical description of the ultrasound system, refer to Appendix A in this manual.

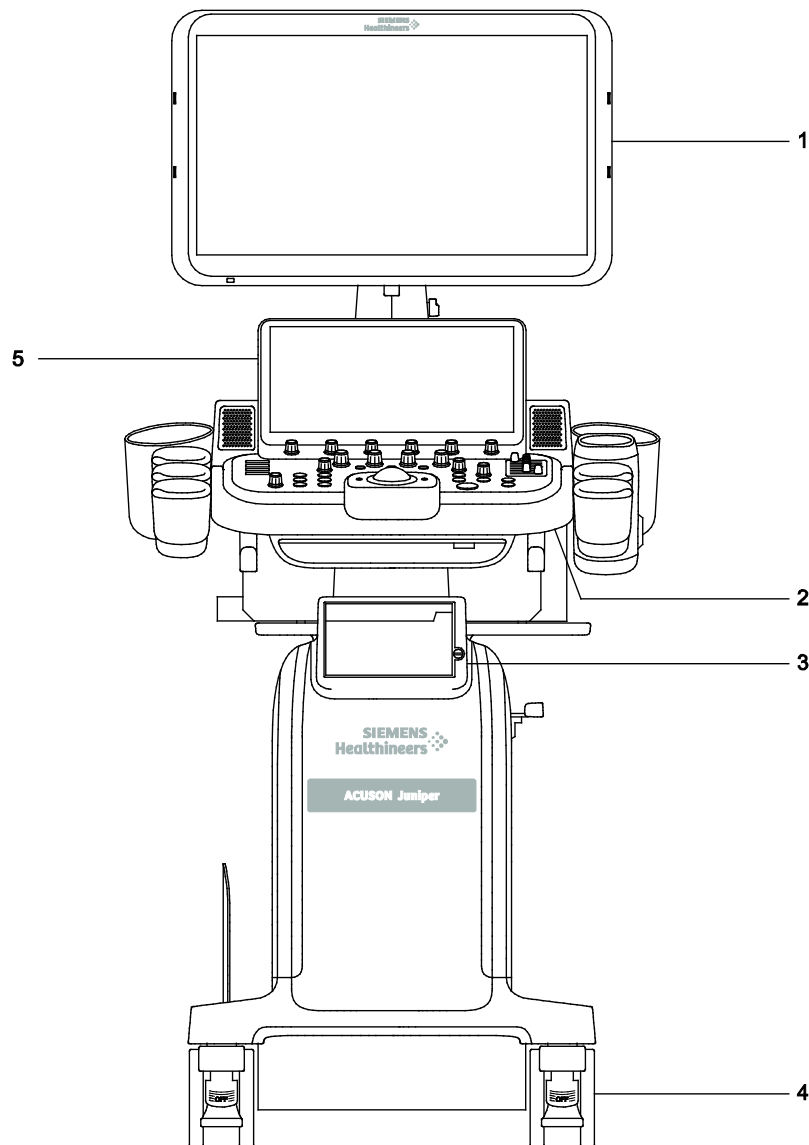


Example of the ultrasound system, front view.

- 1 User-adjustable monitor
- 2 Control panel
- 3 Black and white printer
- 4 Wheel
- 5 Touch screen

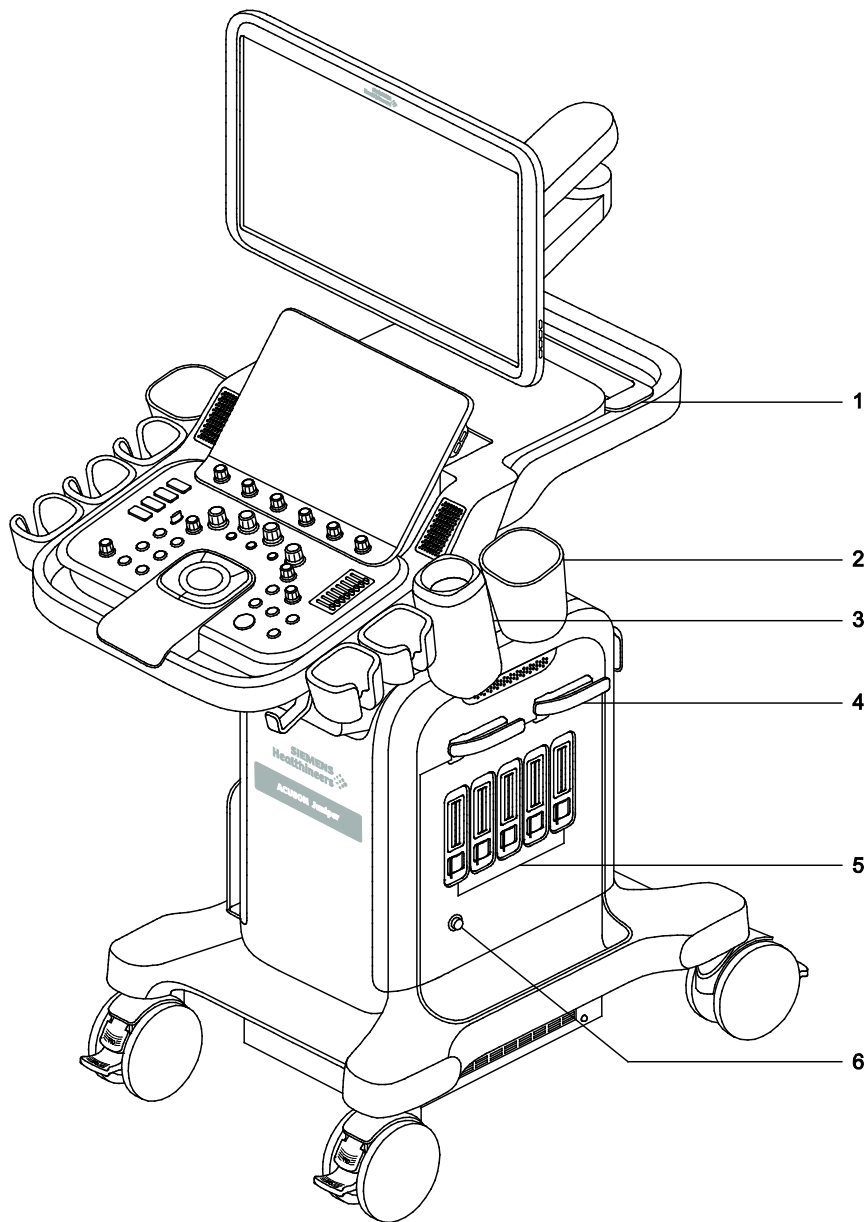


Example of the ultrasound system, right front view.

- 1 Rear basket
- 2 Holder for large bottle of gel
- 3 Gel warmer
- 4 Transducer cable hanger
- 5 Transducer ports
- 6 Transducer port for continuous wave (pencil) transducers

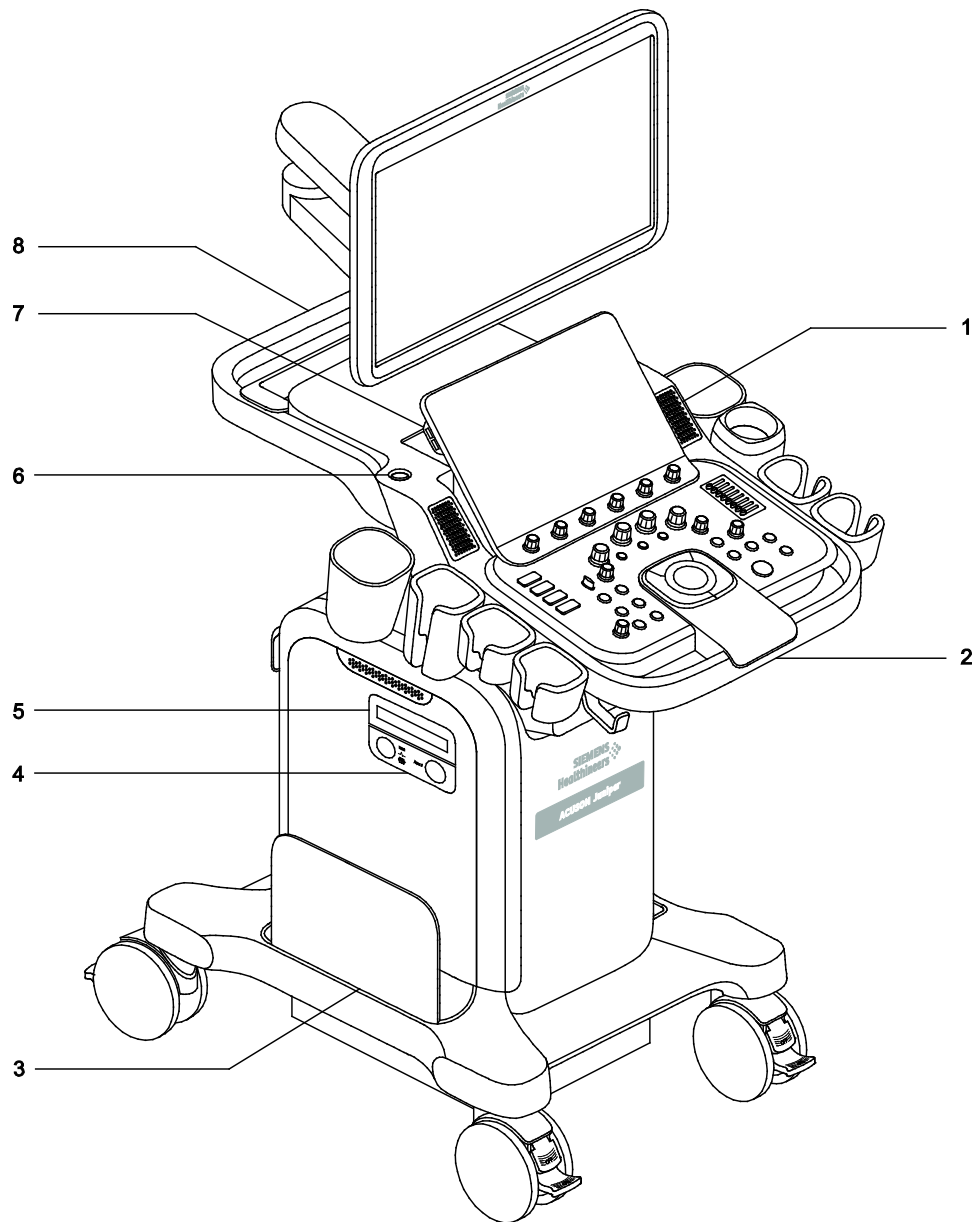

Example of the ultrasound system, left front view.

- 1 Speaker
- 2 Front handle
- 3 Side storage bin
- 4 Electrocardiogram (ECG) module
- 5 CD/DVD read/write drive
- 6 Power on/off
- 7 Touch screen with USB ports
- 8 Rear handle

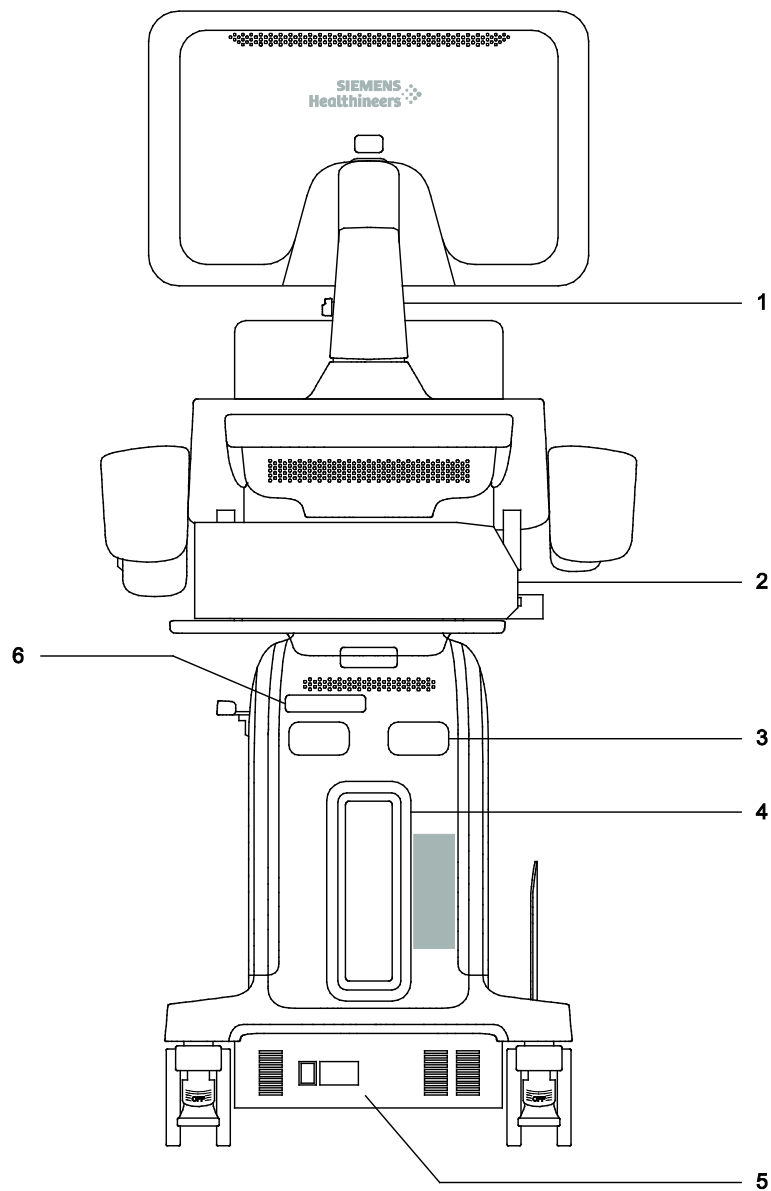

Example of the ultrasound system, back view.

- 1 Adjustable arm
- 2 Printer or DVR
- 3 Power cable hanger
- 4 Input/output panel
- 5 Power supply unit
- 6 Air filter

Intended Use

-  **WARNING:** The analysis of results from an ultrasound examination requires that you are trained in the interpretation of diagnostic ultrasound studies and are qualified to make clinical diagnoses.
-  **WARNING:** To avoid the risk of human injury and permanent damage to the ultrasound system, do not use the ultrasound system, transducers, or ultrasound components and accessories during magnetic resonance imaging (MRI) or in an environment with magnetic resonance (MR) equipment.
-  **WARNING:** The ultrasound system is not intended for ophthalmic use.
-  **Caution:** In the United States of America, federal law restricts this device to sale or use by, or on the order of, a physician.
-  **Caution:** Ultrasound is used as an imaging aid, but may have further restrictions specific to in-vitro fertilization (IVF), chorionic villus sampling (CVS), and percutaneous umbilical cord blood sampling (PUBS) procedures. Observe local laws and regulations.

Indications for Use Statement

Product	Indications for Use Statement
ACUSON Juniper Ultrasound System	The ACUSON Juniper ultrasound imaging system is intended to provide images of, or signals from, inside the body by an appropriately trained healthcare professional in a clinical setting for the following applications: Abdominal, Obstetrics, Gynecology, Small Parts, Pediatric, Neonatal, Vascular, Urology, Echocardiography, Musculoskeletal, and Intraoperative applications using different ultrasound transducers for different applications. The system supports the Ultrasound-Derived Fat Fraction (UDFF) measurement tool to report an index that can be useful as an aid to a physician managing adult patients with hepatic steatosis. The system also provides the ability to measure anatomical structures and provides analysis packages that provide information used by a physician for clinical diagnostic purposes.
ACUSON Juniper Select Ultrasound System	The ACUSON Juniper Select ultrasound imaging system is intended to provide images of, or signals from, inside the body by an appropriately trained healthcare professional in a clinical setting for the following applications: Abdominal, Obstetrics, Gynecology, Small Parts, Pediatric, Vascular, Urology, Echocardiography, and Musculoskeletal applications using different ultrasound transducers for different applications. The system also provides the ability to measure anatomical structures and provides analysis packages that provide information used by a physician for clinical diagnostic purposes.
<i>syngo</i> Arterial Health Package (AHP)*	The Arterial Health Package (AHP) software provides the physician with the capability to measure Intima Media Thickness and the option to reference normative tables that have been validated and published in peer-reviewed studies. The information is intended to provide the physician with an easily understood tool for communicating with patients regarding the state of their cardiovascular system.

***Note:** This feature should be utilized according to the "ASE Consensus Statement; Use of Carotid Ultrasound to Identify Subclinical Vascular Disease and Evaluate Cardiovascular Disease Risk: A Consensus Statement from the American Society of Echocardiography Carotid Intima-Media Thickness Task Force. Endorsed by the Society for Vascular Medicine."

The transducer is intended to provide images of, or signals from, inside the body by an appropriately trained healthcare professional in a clinical setting.

Operating Modes

- 2D-mode
 - 2D-mode with Harmonic imaging
 - 2D-mode with Simultaneous imaging for biplane transducers*
- Color flow Doppler
 - Color (velocity)
 - Power (energy)
- Pulsed Wave Doppler
 - Pulsed Wave Doppler Tissue Imaging
 - High Pulsed Repetition Frequency Pulsed Wave Doppler
- Continuous Wave Doppler
 - Steerable Continuous Wave Doppler for phased array transducers
 - Auxiliary Continuous Wave Doppler for pencil transducers
- M-mode
 - M-mode with Harmonic imaging
 - Anatomical M-mode
- 3D/4D Volume Imaging

Combined Modes

- 2D-mode with Color
- 2D-mode with Doppler
- 2D-mode with Color and Doppler
- 2D-mode with M-mode
- 2D-mode with M-mode and Color
- 2D-mode with Elastography*
- 2D-mode with Contrast Agent Imaging
- 2D-mode with Anatomical M-mode

*Available only on the ACUSON Juniper system.

Patient Target Groups

- Adult
- Pediatric
- Fetal

Intended Use Environment

The ultrasound system is a mobile device intended for use in clinical environments.

See also: For information about the environmental requirements, refer to Appendix A in this manual.

Intended Users

User	Interaction with Ultrasound Equipment	Expected Experience and Other Characteristics
Sonographer	- Acquires diagnostic views of anatomy, blood flow, and related pathology - Performs measurements and analysis of the acquired images - Prepares exam data for review and interpretation by a qualified physician	- Ranges from novices (for example, students) to advanced practitioners with certification in multiple subspecialties - Educated in anatomy, physiology, patient care, and identification of pathology in ultrasound images - Many sonographers have a Bachelor's degree; some have advanced degrees in related health care subjects
Physician	- Performs ultrasound exams - Interprets exam data - Writes and assembles exam findings in a report	- Medical doctor - Trained in ultrasound imaging techniques - Skilled in interpreting ultrasound exam data
Radiologist	- Performs ultrasound exams - Performs image-guided access and biopsy procedures - Interprets exam data - Writes and assembles exam findings in a report	- Medical doctor - Expert in diagnostic imaging, including computer tomography, magnetic resonance imaging, X-ray, ultrasound, and nuclear medicine - Advanced training in imaging physics with typically two to six years of post-doctoral training in the field of radiology; some specialize in diagnostic ultrasound
Cardiologist	- Performs ultrasound exams - Interprets exam data - Writes and assembles exam findings in a report	- Medical doctor - Expert in diagnostic imaging, including computer tomography, magnetic resonance imaging, X-ray, ultrasound, and nuclear medicine - Advanced training in imaging physics with typically two to six years of post-doctoral training in the field of cardiology; some specialize in diagnostic ultrasound
System Administrator and Customer Service Engineer	Configures the ultrasound system for use in a networked environment	- A system administrator is an individual within your organization who is designated to set up system parameters to connect the ultrasound system or workstation to a picture archiving and communication system (PACS). - A customer service engineer is an individual representing the manufacturer of the ultrasound system who configures the ultrasound system during initial installation, supports troubleshooting activities, and repairs the system.

Contraindications

There are no known contraindications for diagnostic ultrasound.

Transducers and Intended Applications

⚠ WARNING: The ultrasound system is not intended to be used with high-frequency surgical equipment. Use of transducers in conjunction with high-frequency surgical equipment may result in an increased risk of harm to the patient in the event of a defect in the high-frequency surgical neutral electrode connection.

Only the following transducers from Siemens Healthineers are compatible with your ultrasound system.


Curved and Linear Array Transducers



The following transducers are available on the ACUSON Juniper ultrasound system:



- 5C1, 6C1, 7C2, 12L3, 11L4, 13L4, 14L4, 16L4, 18H5, 11M3, L10-5v



The following transducers are available on the ACUSON Juniper Select ultrasound system:




- 5C1, 12L3, 11L4, L10-5v

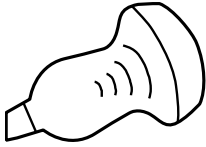Transducer Name	Modes of Operation ¹	Intended Applications	Exam Types on the System
5C1 	2D, C, D, M	▪ Abdominal ▪ Abdominal Difficult ▪ Renal ▪ Bowel ▪ Focused Assessment with Sonography for Trauma ▪ Lung ▪ Obstetrics ▪ Early Obstetrics ▪ Advanced Obstetrics ▪ Fetal Echocardiography ▪ Gynecology ▪ Peripheral Venous ▪ Peripheral Arterials ▪ Pelvis ▪ Spine	▪ Abd ▪ Abd-Difficult ▪ Renal ▪ Bowel ▪ FAST ▪ Lung ▪ OB ▪ Early OB ▪ OB(Adv) ▪ Fetal Echo ▪ GYN ▪ Venous ▪ Arterial ▪ Pelvis ▪ Spine

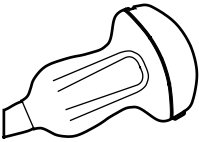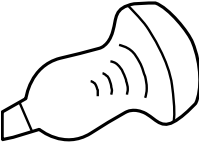Transducer Name	Modes of Operation ¹	Intended Applications	Exam Types on the System
6C1 	2D, C, D, M	▪ Abdominal ▪ Abdominal Difficult ▪ Renal ▪ Bowel ▪ Focused Assessment with Sonography for Trauma ▪ Lung ▪ Obstetrics ▪ Early Obstetrics ▪ Advanced Obstetrics ▪ Fetal Echocardiography ▪ Gynecology ▪ Peripheral Venous ▪ Peripheral Arterials ▪ Pelvis ▪ Spine	▪ Abd ▪ Abd-Difficult ▪ Renal ▪ Bowel ▪ FAST ▪ Lung ▪ OB ▪ Early OB ▪ OB(Adv) ▪ Fetal Echo ▪ GYN ▪ Venous ▪ Arterial ▪ Pelvis ▪ Spine
7C2 	2D, C, D, M	▪ Abdominal ▪ Renal ▪ Bowel ▪ Focused Assessment with Sonography for Trauma ▪ Lung ▪ Obstetrics ▪ Early Obstetrics ▪ Advanced Obstetrics ▪ Fetal Echocardiography ▪ Gynecology ▪ Spine ▪ Pediatric Abdomen ▪ Pelvis	▪ Abd ▪ Renal ▪ Bowel ▪ FAST ▪ Lung ▪ OB ▪ Early OB ▪ OB(Adv) ▪ Fetal Echo ▪ GYN ▪ Spine ▪ Ped Abd ▪ Pelvis

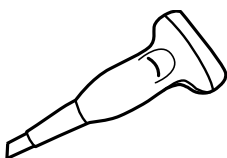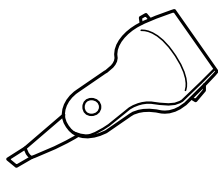Transducer Name	Modes of Operation ¹	Intended Applications	Exam Types on the System
12L3 	2D, C, D, M	▪ Bowel ▪ Focused Assessment with Sonography for Trauma ▪ Lung ▪ Obstetrics ▪ Breast ▪ Testis ▪ Penile ▪ Thyroid ▪ Pediatric Hip Joint ▪ Carotid ▪ Peripheral Venous ▪ Peripheral Arterials ▪ Spine ▪ Musculoskeletal ▪ Digital ▪ Nerve	▪ Bowel ▪ FAST ▪ Lung ▪ OB ▪ Breast ▪ Testis ▪ Penile ▪ Thyroid ▪ Ped Hip ▪ Carotid ▪ Venous ▪ Arterial ▪ Spine ▪ MSK ▪ Digital ▪ Nerve
11L4 	2D, C, D, M	▪ Bowel ▪ Focused Assessment with Sonography for Trauma ▪ Lung ▪ Obstetrics ▪ Breast ▪ Testis ▪ Penile ▪ Thyroid ▪ Pediatric Hip Joint ▪ Carotid ▪ Peripheral Venous ▪ Peripheral Arterials ▪ Spine ▪ Musculoskeletal ▪ Digital ▪ Nerve	▪ Bowel ▪ FAST ▪ Lung ▪ OB ▪ Breast ▪ Testis ▪ Penile ▪ Thyroid ▪ Ped Hip ▪ Carotid ▪ Venous ▪ Arterial ▪ Spine ▪ MSK ▪ Digital ▪ Nerve

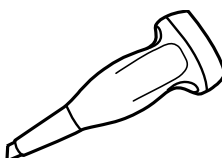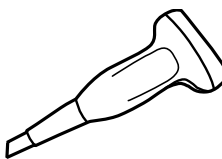Transducer Name	Modes of Operation ¹	Intended Applications	Exam Types on the System
13L4 	2D, C, D, M	▪ Bowel ▪ Focused Assessment with Sonography for Trauma ▪ Lung ▪ Obstetrics ▪ Breast ▪ Testis ▪ Penile ▪ Thyroid ▪ Pediatric Hip Joint ▪ Carotid ▪ Peripheral Venous ▪ Peripheral Arterials ▪ Spine ▪ Musculoskeletal ▪ Digital ▪ Nerve	▪ Bowel ▪ FAST ▪ Lung ▪ OB ▪ Breast ▪ Testis ▪ Penile ▪ Thyroid ▪ Ped Hip ▪ Carotid ▪ Venous ▪ Arterial ▪ Spine ▪ MSK ▪ Digital ▪ Nerve
14L4 	2D, C, D, M	▪ Bowel ▪ Focused Assessment with Sonography for Trauma ▪ Lung ▪ Obstetrics ▪ Breast ▪ Testis ▪ Penile ▪ Thyroid ▪ Pediatric Hip Joint ▪ Carotid ▪ Peripheral Venous ▪ Peripheral Arterials ▪ Spine ▪ Musculoskeletal ▪ Digital ▪ Nerve	▪ Bowel ▪ FAST ▪ Lung ▪ OB ▪ Breast ▪ Testis ▪ Penile ▪ Thyroid ▪ Ped Hip ▪ Carotid ▪ Venous ▪ Arterial ▪ Spine ▪ MSK ▪ Digital ▪ Nerve

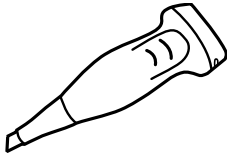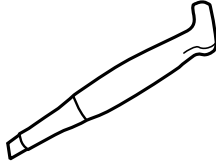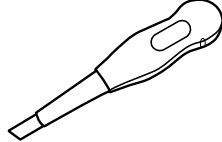Transducer Name	Modes of Operation ¹	Intended Applications	Exam Types on the System
16L4 	2D, C, D, M	▪ Bowel ▪ Focused Assessment with Sonography for Trauma ▪ Lung ▪ Breast ▪ Testis ▪ Penile ▪ Thyroid ▪ Pediatric Hip Joint ▪ Carotid ▪ Musculoskeletal ▪ Digital ▪ Nerve	▪ Bowel ▪ FAST ▪ Lung ▪ Breast ▪ Testis ▪ Penile ▪ Thyroid ▪ Ped Hip ▪ Carotid ▪ MSK ▪ Digital ▪ Nerve
18H5 	2D, C, D, M	▪ Testis ▪ Penile ▪ Thyroid ▪ Musculoskeletal ▪ Digital ▪ Nerve ▪ Intraoperative Abdomen ▪ Carotid	▪ Testis ▪ Penile ▪ Thyroid ▪ MSK ▪ Digital ▪ Nerve ▪ IO Abd ▪ Carotid
11M3 	2D, C, D, M	▪ Pediatric Abdomen ▪ Neonatal Echocardiography ▪ Neonatal Head ▪ Carotid ▪ Pediatric Echocardiography	▪ Ped Abd ▪ Neo Echo ▪ Neo Head ▪ Carotid ▪ Ped Echo






Phased Array Transducers

The following transducers are available on the ACUSON Juniper ultrasound system:

- 5P1, 8V4, 10V4, 5VT, P4-2

The following transducers are available on the ACUSON Juniper Select ultrasound system:

- 5P1, P4-2

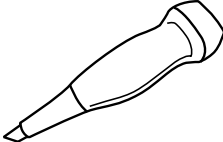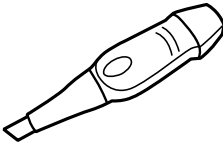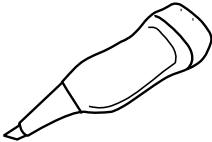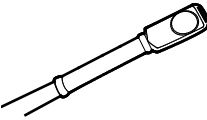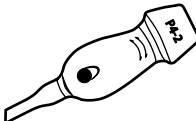Transducer Name	Modes of Operation ¹	Intended Applications	Exam Types on the System
5P1 	2D, C, D, M, CW	▪ Abdominal ▪ Renal ▪ Focused Assessment with Sonography for Trauma ▪ Lung ▪ Transcranial Doppler ▪ Cardiac ▪ Pediatric Echocardiography	▪ Abd ▪ Renal ▪ FAST ▪ Lung ▪ TCD ▪ Cardiac ▪ Ped Echo
8V4 	2D, C, D, M, CW	▪ Pediatric Abdomen ▪ Neonatal Echocardiography ▪ Neonatal Head ▪ Pediatric Echocardiography	▪ Ped Abd ▪ Neo Echo ▪ Neo Head ▪ Ped Echo
10V4 	2D, C, D, M, CW	▪ Pediatric Abdomen ▪ Neonatal Echocardiography ▪ Neonatal Head ▪ Pediatric Echocardiography	▪ Ped Abd ▪ Neo Echo ▪ Neo Head ▪ Ped Echo
5VT 	2D, C, D, M, CW	▪ Transesophageal Echocardiography	▪ TEE
P4-2 	2D, C, D, M, CW	▪ Abdominal ▪ Transcranial Doppler ▪ Cardiac	▪ Abd ▪ TCD ▪ Cardiac




Endocavity Transducers

The following transducers are available on the ACUSON Juniper ultrasound system:

- 9MC3, 10MC3, 10BCC3

The following transducers are available on the ACUSON Juniper Select ultrasound system:

- 10MC3

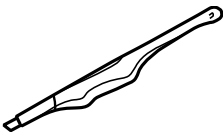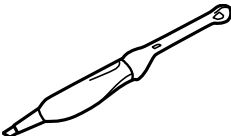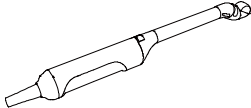Transducer Name	Modes of Operation ¹	Intended Applications	Exam Types on the System
9MC3 	2D, C, D, M	▪ Obstetrics ▪ Early Obstetrics ▪ Advanced Obstetrics ▪ Gynecology ▪ Prostate	▪ OB ▪ Early OB ▪ OB(Adv) ▪ GYN ▪ Prostate
10MC3 	2D, C, D, M	▪ Obstetrics ▪ Early Obstetrics ▪ Advanced Obstetrics ▪ Gynecology ▪ Prostate	▪ OB ▪ Early OB ▪ OB(Adv) ▪ GYN ▪ Prostate
10BCC3 	2D, C, D, M	▪ Prostate	▪ Prostate


fourSight 4D Transducers

The following transducers are available on the ACUSON Juniper ultrasound system:

- 9VC2, 8VC3, 9VE3, 9VE4

The following transducers are available on the ACUSON Juniper Select ultrasound system:

- 9VC2, 8VC3

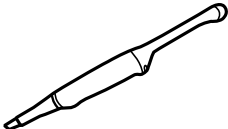Transducer Name	Modes of Operation ¹	Intended Applications	Exam Types on the System
8VC3 	2D, C, D, M, Volume Imaging	▪ Abdominal ▪ Obstetrics ▪ Early Obstetrics ▪ Advanced Obstetrics ▪ Gynecology ▪ Fetal Echocardiography ▪ Pelvic Floor	▪ Abd ▪ OB ▪ Early OB ▪ OB(Adv) ▪ GYN ▪ Fetal Echo ▪ Pelvic Floor
9VC2 	2D, C, D, M, Volume Imaging	▪ Abdominal ▪ Obstetrics ▪ Early Obstetrics ▪ Advanced Obstetrics ▪ Fetal Echocardiography ▪ Gynecology ▪ Pelvic Floor	▪ Abd ▪ OB ▪ Early OB ▪ OB(Adv) ▪ Fetal Echo ▪ GYN ▪ Pelvic Floor
9VE3 	2D, C, D, M, Volume Imaging	▪ Obstetrics ▪ Early Obstetrics ▪ Advanced Obstetrics ▪ Gynecology ▪ Pelvic Floor	▪ OB ▪ Early OB ▪ OB(Adv) ▪ GYN ▪ Pelvic Floor
9VE4 	2D, C, D, M, Volume Imaging	▪ Obstetrics ▪ Early Obstetrics ▪ Advanced Obstetrics ▪ Gynecology ▪ Pelvic Floor	▪ OB ▪ Early OB ▪ OB(Adv) ▪ GYN ▪ Pelvic Floor



CW Transducers

The following transducers are available on the ACUSON Juniper ultrasound system:

- CW2, CW5, CW8

The following transducers are available on the ACUSON Juniper Select ultrasound system:

- CW2, CW5

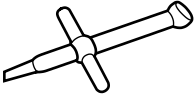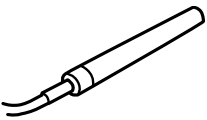Transducer Name	Modes of Operation ¹	Intended Applications	Exam Types on the System
CW2 	CW	▪ Cardiac ▪ Pediatric Echocardiography	▪ Cardiac ▪ Ped Echo
CW5 	CW	▪ Carotid ▪ Transcranial Doppler ▪ Pediatric Echocardiography	▪ Carotid ▪ TCD ▪ Ped Echo
CW8 	CW	▪ Carotid ▪ Peripheral Arterials	▪ Carotid ▪ Arterial

Legend

¹ Modes of Operation	Includes one or more of the following system operating modes:
2D (brightness mode)	2D-mode, 2D-mode with Harmonics imaging, 2D-mode with Simultaneous imaging (for biplane transducers), 2D-mode with Contrast Agent Imaging
C (color flow Doppler)	Color flow Doppler, Color (velocity), Power (energy), 2D-mode with Color
D (Doppler)	Pulsed Wave Doppler, Pulsed Wave Doppler Tissue Imaging, High Pulsed Repetition Frequency Pulsed Wave Doppler, 2D-mode with Doppler, 2D-mode with Color and Doppler
M (motion mode)	M-mode, M-mode with Harmonics imaging, 2D-mode with M-mode, 2D-mode with M-mode and Color, Anatomical M-mode, 2D-mode with Anatomical M-mode
CW (continuous wave Doppler)	Steerable Continuous Wave Doppler (for phased array transducers), Auxiliary Continuous Wave Doppler (for continuous wave [pencil] transducers)
Volume Imaging	3D volume imaging, 4D volume imaging

Principles of Operation

Diagnostic ultrasound uses high-frequency (above the audible range) sound waves to produce an image of anatomical structures within the body. Electrical pulses vibrate ceramics within a transducer to transmit the sound waves into the body. Sound waves travel through body tissue at approximately 1,540 meters per second and reflect back as echoes to the transducer at each point of change in tissue density, for example, at the border of two organs in the body. These return signals provide information about the acoustic properties of the tissue. The time to receive the echo in microseconds indicates a distance into the body. Structures furthest from the transducer surface require more time to return a signal than structures closer to the transducer surface. The strength and position of each signal indicates a point of varying intensity (brightness). The distances and intensities are processed and displayed on a screen to form a two-dimensional (2D) image.

Clinical Benefits

Ultrasound has a positive impact on the health of an individual as a non-ionizing diagnostic imaging tool. The device provides real-time anatomic and kinetic visualization of internal organs, tissue structures, blood flow, and Doppler-enabled hemodynamics. Ultrasound is the accepted imaging modality for the safe visualization of fetal anatomy and monitoring of fetal development.

The mobility and cost effectiveness of ultrasound increases access to patients in clinical environments. Patient benefits also include a safer examination compared to imaging modalities using ionizing radiation.

Real-time ultrasound imaging provides guidance during interventional procedures, for example, needle biopsies and fluid aspirations; in the absence of ultrasound, these procedures are conducted without visualization or with ionizing radiation.


Image Screen Layout

The monitor on the ultrasound system displays clinical images together with important operating parameters and patient data.

Screen Saver

The screen saver feature freezes the ultrasound system and replaces the active display with a screen saver after the system has been inactive for a specified number of minutes. A user login is required to exit the screen saver and access the system.

Note: The screen saver feature is not available when the system is in external video playback or during the Biopsy function.

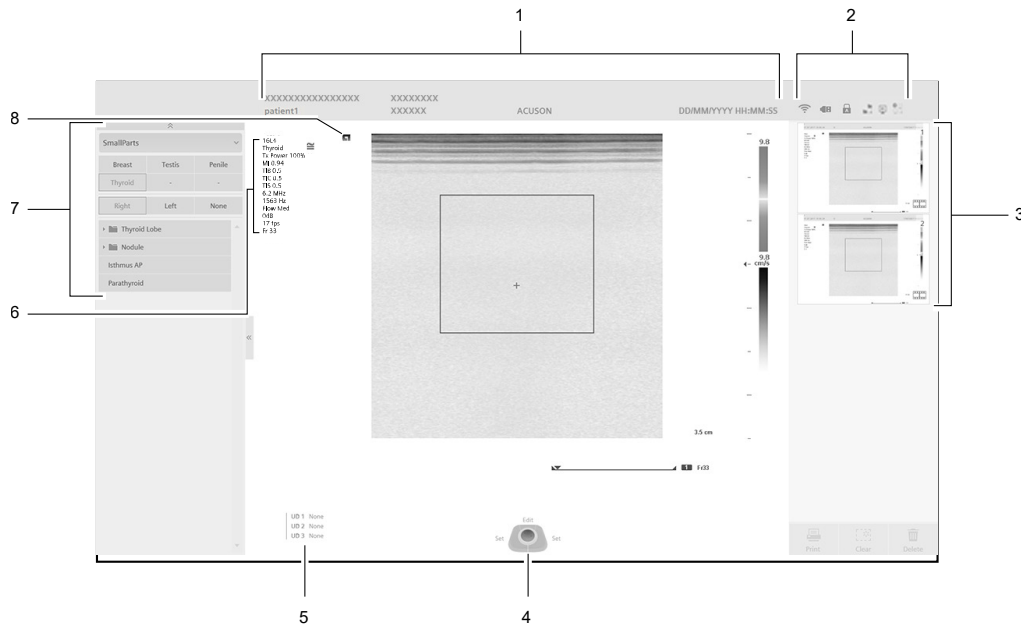

Example of a typical image screen.

- 1 Patient banner
Displays patient information, operator, institution, the system date and time.
- 2 System status indicators
Displays information about the ultrasound system, for example, the signal strength of a wireless network and the status of printing, video recording, and transferring data.
- 3 Thumbnail panel
Displays reference images of stored clips, images, and volumes.
- 4 Trackball legend
Identifies the functions currently assigned to the trackball and surrounding controls on the control panel. In this manual, the trackball keys are named **SET**.

- 5 User-defined key legend
Identifies the function currently assigned to user-defined keys.
- 6 Imaging parameters
Displays a list of imaging settings for the active operating mode.
- 7 Left panel
Provides functions for workflow protocols, reports, and for specialty features, for example, Point Shear Wave Elastography.
During the measurement function, displays measurement labels and modifiers.
- 8 Transducer orientation indicator
Indicates the direction of the scan.

System Status Indicators


The system status is always displayed to provide information about the system.

Icon	Description
	Indicates a software download to the ultrasound system is in progress.
	Indicates remote service is in progress.
	Indicates the system is connected to a wireless network and indicates the signal strength.
	Indicates a wireless network is not found.
	Indicates a wireless network error.
	Indicates the system is connected to a wired network.
	Indicates the system is trying to connect to a wired network.
	Indicates an error occurred while connecting to a wired network.

Icon	Description
	Indicates printing is in progress.
	Indicates an error occurred during printing.
	Indicates a data transfer is in progress.
	Indicates a data transfer error.
	Indicates the status of video recording.
	Indicates the status of video paused.
	Indicates a USB device is connected to the system.
	Indicates the CAPS LOCK key is activated.
	Indicates the system is connected to the HDMI.

About the Control Panel

The control panel provides access to system features and workflow tasks. The keys and controls are logically arranged to require a minimum of hand and eye movements. The control panel also includes customizable keys.

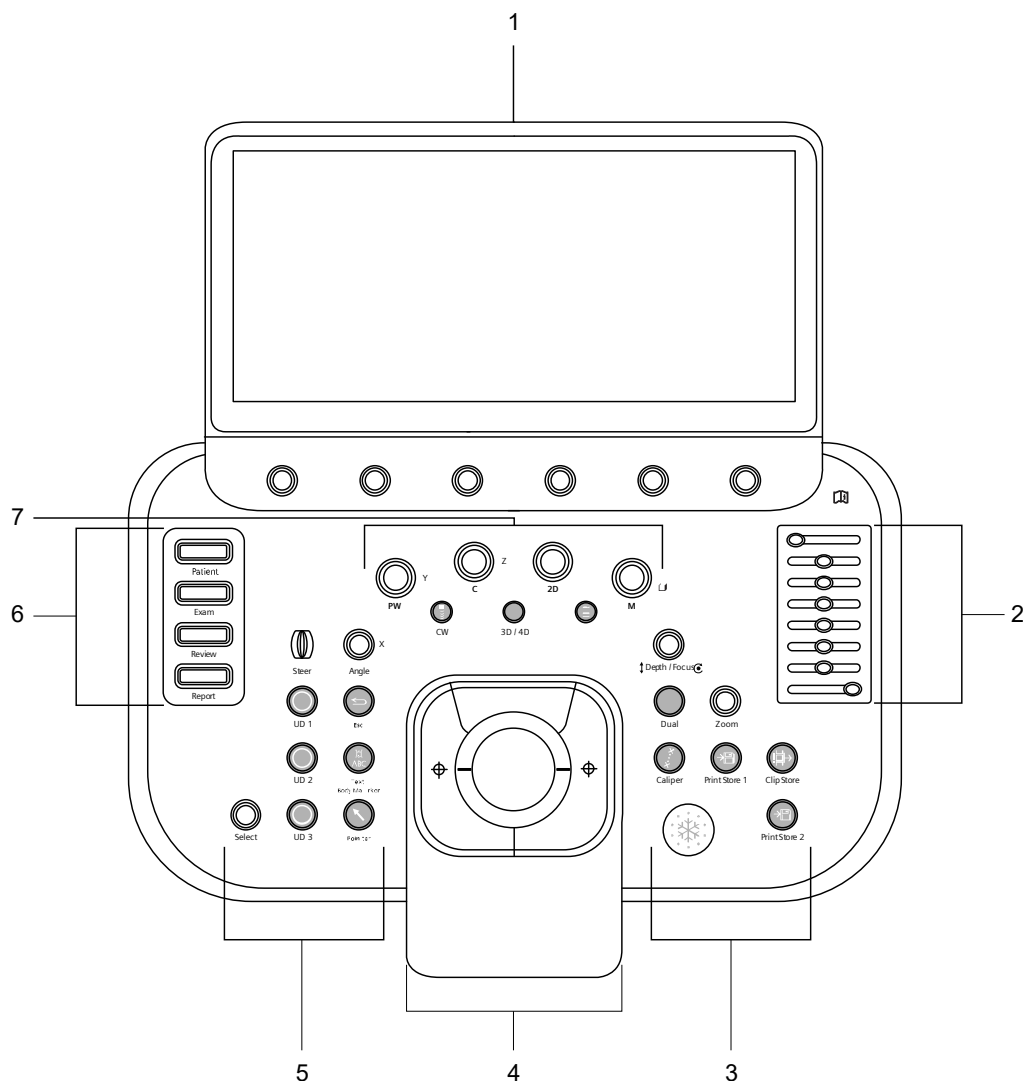

Example of control panel on the ultrasound system.

- 1 Touch screen with touch screen controls
Use the six controls to change imaging parameter settings, select a measurement method, or perform review functions. The screen displays the functions assigned to each key.
- 2 Depth Gain Compensation controls
Slide the controls left or right to manually adjust the depth-gain compensation.

3 Image acquisition and documentation controls

Use these controls to freeze the image, adjust the depth or zoom factor, perform measurements, and save or print an image or clip.

4 Trackball controls: the trackball, **SET** keys and **UPDATE**

Use the trackball to position or resize objects on the image screen, for example, the 2D FOV (field of view), the color ROI (region of interest), a Doppler cursor, and measurement markers.

The function of all trackball controls depends on the active task. Use the **SET** keys to anchor the measurement markers, and use the **UPDATE** key to activate or deactivate an operating mode, acquire imaging volume or clips, or to refresh a 2D-mode image.

5 Doppler optimization, annotation controls, and user-defined controls

Use the Doppler optimization controls to adjust the color or correct the Doppler flow angle.

Press the annotation controls to activate the annotation function including text, body markers, and pointer.

Press the user-defined controls to activate the function specified in the configuration settings.

6 Workflow controls

Use these controls to activate workflow tasks, for example, registering a patient, displaying a report.

7 Operating mode controls

Press the control to activate the required operating mode. Rotate the control to adjust the gain for the active mode.

Control Panel Lighting

When the ultrasound system is in use, controls and keys on the control panel are supported by backlighting. The brightness of the backlighting increases when a control is active.


Power Control



Press this control to power on or power off the ultrasound system.

Workflow Controls


Note: Press the workflow control to activate or exit the function.

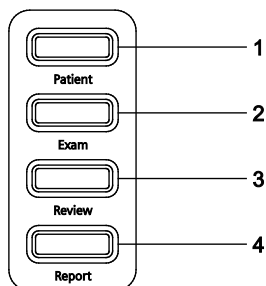

Example of workflow controls next to the alphanumeric keyboard.

- 1 **PATIENT** displays the patient registration form. You can enter patient data and select an exam type and transducer.
- 2 **EXAM** displays a list on the touch screen of transducers connected to the ultrasound system and the supported exam types for the selected transducer.
- 3 **REVIEW** displays the review screen for the current exam. You can view, print, export, or delete saved images and clips.
- 4 **REPORT** displays the report for the current exam. You can review and edit results, add images to the report, preview and print the report, and export the report in PDF format.


Operating Mode Controls

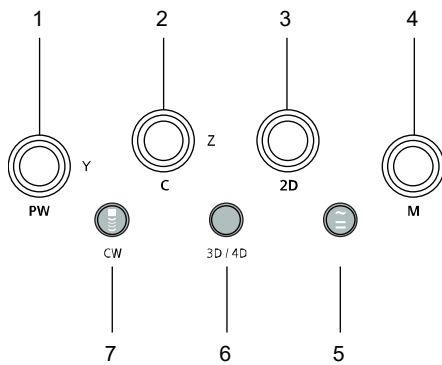

Example of operating mode controls.

- 1 **PW**
Press this control to activate pulsed wave (PW) Doppler and display the Doppler spectrum and 2D image as specified in the configuration settings.
Rotate this control to adjust the Doppler gain. Use the touch screen selections to activate Doppler Tissue Imaging (DTI).
- 2 **C**
Press this control to activate color Doppler.
Rotate this control to adjust the color gain. Use the touch screen selections to activate additional Doppler options, for example, color or power
- 3 **2D**
Press this control to activate 2D-mode or to exit another operating mode.
Rotate this control to adjust the 2D gain.
- 4 **M**
Press this control to activate M-mode and display the M-mode sweep and 2D image as specified in the configuration settings.
Rotate this control to adjust the M-mode gain.
- 5 **AUTO TEQ**
Press this control to optimize the uniformity of image brightness in the field of view.
- 6 **3D / 4D**
Press this control to activate volume imaging.
- 7 **CW**
Press this control to activate continuous wave (CW) Doppler or steerable continuous wave (SCW) Doppler function for phased array transducers.

Image Acquisition and Documentation Controls

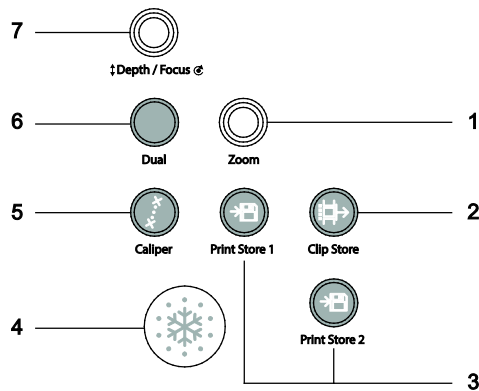

Example of image acquisition and documentation controls.

1 ZOOM

Rotate this control to change the image magnification.
Press the control to activate a high-density zoom image on the screen.

2 CLIP STORE

Press this control to acquire and save a clip or static image to the system's hard disk as specified in the configuration settings. The thumbnail panel displays the saved data.

During Stress Echo, this control saves a clip as specified in the configuration setting for a continuous protocol or alternates beginning and ending acquisition for a non-continuous protocol.

3 PRINT STORE 1, PRINT STORE 2

Press this control to acquire a static image. The system saves the image to the system's hard disk, sends the image to a printer, or both, as specified in the configuration settings.

4 FREEZE

Press this control to freeze an image, sweep, or spectrum or restore real-time imaging.

During system freeze, you can review or save an image or clip, perform measurements, or add annotations or body markers.

5 CALIPER

Press this control to activate or deactivate the measurement function.


6 DUAL

Press this control to activate dual imaging.

7 DEPTH / FOCUS

Rotate or push this control to adjust the imaging depth as specified in the configuration settings. Press the control to add multiple transmit focal zones. Rotate or push the control to reposition the transmit focal zones in the image as specified in the configuration settings.

Optimization Controls

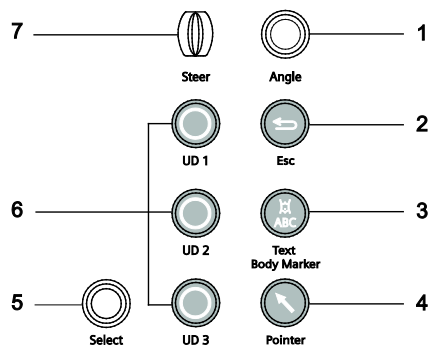

Example of optimization controls.

1 ANGLE

During Doppler, press this control to adjust the flow angle in 60° increments (60°, 0°, and -60°). Rotate the control to adjust the flow angle in 1° increments.

During Panoramic imaging, rotate this control to rotate a composite image.

During annotation with a body marker, press this control to rotate the transducer marker in 90° increments. Rotate the control to adjust the angle of the transducer marker in 15° increments.

2 ESC

Exits the currently displayed mode, function, or page and reactivates the previous mode, function, or page.

3 TEXT BODY MARKER

Activates the text, body marker, and arrows on the image.

4 POINTER

Press the pointer control to access the controls or thumbnails on the image screen.

5 SELECT

Press or rotate this control to make selections from the screen.

6 UD 1 / UD 2 / UD 3

Activates the function assigned to this control specified in the configuration settings.

7 STEER


Positions the M-mode or Doppler cursor or the color region of interest.

For linear array transducers, steers the 2D image. The amount of steering is dependent on the type of transducer.

Retractable Keyboard

⚠ WARNING: Carefully extend or retract the keyboard. A potential pinch point between the control panel and keyboard can create a risk for injury to the user.

The retractable alphanumeric keyboard attaches below the control panel. The retractable keyboard is an ultrasound system option, not a standard feature.

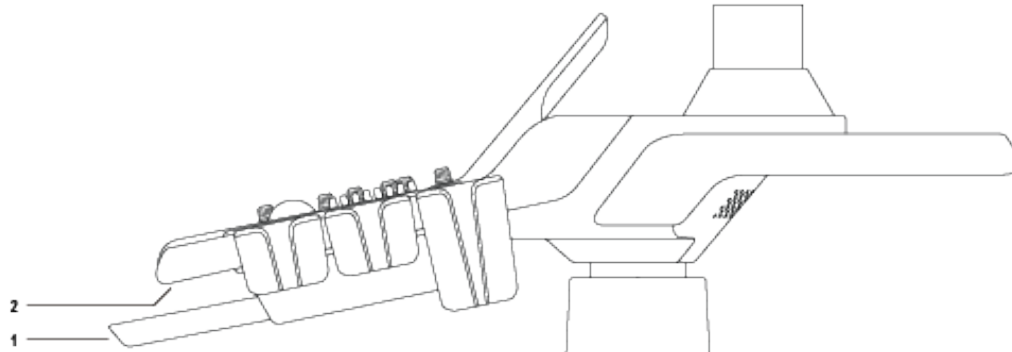

Example of the retractable keyboard on the ultrasound system.

- 1 Push point
- 2 Bottom of the control panel

To access the keyboard:

- To extend the keyboard, press the push point until you hear an audible click.
- The keyboard releases from the control panel and extends for use.
- To retract the keyboard, push the keyboard toward the ultrasound system until you hear an audible click.

Special Keys

Key Name	Description
Text A, Text B	Inserts text on the image screen at the cursor location.
Arrow	Inserts an arrow on the image screen at the cursor location.
END	Closes the study as specified in the configuration settings.

Special Characters

To access special characters on the retractable keyboard:

1. Press and hold the **Alt** key on the keyboard.
2. Press one of the following keys on the keyboard:

Spanish		German		French	
Key	Character	Key	Character	Key	Character
1	í	p	Ä	n	Ç
/	¿	-	ä	j	ç
i	Ñ	0 (zero)	ö	z	à
9	ñ	[	Ü	x	è
a	á	=	ü	s	é
s	é	'	ß	b	ù
d	í				
f	ó				
g	ú				
Italian		Scandinavian		Other	
Key	Character	Key	Character	Key	Character
z	à	m	Æ	u	Œ
x	è	k	æ	8	œ
s	é	. (period)	Ø	q	α
c	ì	;	ø	w	β
v	ò	,	Å	e	μ
b	ù	l	å	r	θ


Touch Screen Controls

The touch screen has selections for use during an exam for each operating mode and during review.

The following colors indicate the availability of a selection for use.

- Blue indicates the selection is active.
- Gray indicates the selection is available, but not active.
- A dimmed selection is not available for use.

Use the configuration settings to change and rearrange selections on the touch screen.

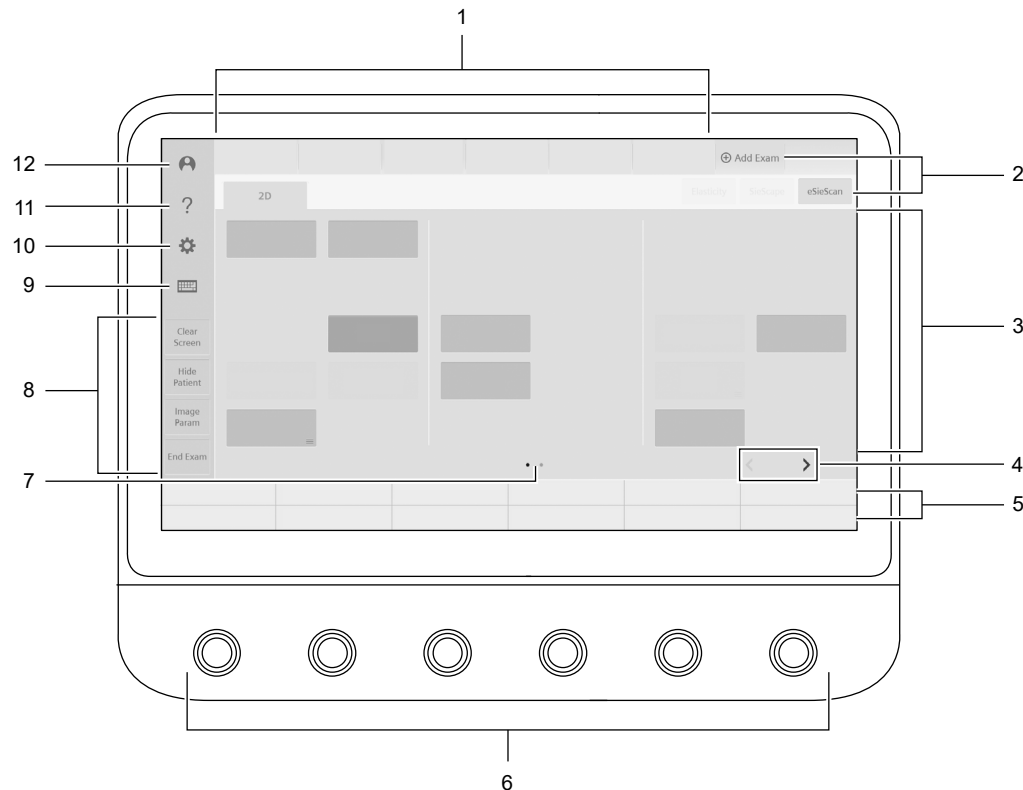

Example of the touch screen layout.

- 1 Selections for the exam type and for transducers connected to the ultrasound system
- 2 Tabs or selections for operating modes, features, or functions
- 3 Selections for the active mode, feature, or function
Tap a selection on the touch screen to activate a feature or function, or to change a setting.
- 4 Additional pages of selections
- 5 Selections for the controls directly below the selections
- 6 Controls corresponding to the touch screen selections directly above the controls
Press, rotate, or press and then rotate a control to make a selection or change a setting
- 7 Number of pages of selections and the active page. The active page is indicated in white.
- 8 Shortcut selections for accessing, hiding, or showing the patient information or image parameters, removing all annotations and pictograms from the screen, and ending the current study
- 9 Selections for accessing the virtual keyboard
- 10 Selections for configuring the ultrasound system
- 11 Information about how to use the ultrasound system
- 12 Options for logging in to the ultrasound system to restrict access to the ultrasound system

Touch Screen Symbols

Symbol	Description
	Indicates a port is available for connecting an array transducer to the ultrasound system.
	Indicates a port is available for connecting a continuous wave (pencil) transducer to the system.
	Accesses additional pages of selections or selections for additional features.
	
	Accesses additional selections.
	Indicates a selection assigned to a control requiring a rotate action.
	Indicates a selection assigned to a control requiring a press action.
	Indicates a package of measurement labels.
	Indicates a group of measurement labels requiring completion before the ultrasound system can calculate a result.
	Indicates a measured or entered value is assigned to the measurement label. During protocols or stress echo, indicates a view is complete.
	Indicates additional options available for the selection.
	- Increases a value. - During stress echo, indicates a protocol or view is temporarily paused. - Enlarges the size of the 2D image.
	
	Temporarily hides the virtual keyboard.
	Rotates the image clockwise for the transesophageal transducer in 90-degree increments.
	Rotates the image counterclockwise for the transesophageal transducer in 90-degree increments.

Touch Screen Keyboard

The touch screen includes an alphanumeric keyboard for entering and editing information in the patient registration form, a report, or on the image.

The on-screen keyboard displays only alphanumeric keys and special function keys related to the active feature. For example, during text annotation, keys specific to text annotation are on the keyboard.

Special Keys and Special Characters

Use the touch screen keyboard to enter text. The keyboard also includes special function keys.

Key Name	Description
English	Accesses special characters specific to a language. Press and hold the key and then select Special .
INS	Inserts text when the text cursor is active.

System Configuration

See also: For detailed information about configuring the ultrasound system, refer to the System Reference manual.

You can customize the system to accommodate your workflow and preferences, for example:

- Language, date and time format, measurement units for height and weight, and audio feedback
- Format of the patient registration form (short or long form)
- Functions assigned to keys on the control panel
- Screen saver, patient banner, and Doppler and M-mode cursor display
- Freeze response, display formats, and other settings by exam type
- Customizable functions assigned to touch screen selections and left panel controls
- Requested and scheduled procedures
- User-defined exam configuration
- Customizable measurement and report settings by exam type
- Annotations and body markers available with each exam type
- Wireless and wired network configuration
- Clip capture type (prospective or retrospective), compression, clip playback, and slide show settings
- Archival settings
- Printer and barcode configuration
- Contrast configuration
- Protocols configuration
- Stress echo configuration
- Security configuration

2 Safety and Care

Operating Safety and Environment.....	2
Serious Incident Reporting	2
System Symbols.....	2
Labels	9
Biohazard Considerations.....	11
Note on Fetal Examinations	12
Acoustic Output	13
Transmit Power Control.....	14
Transmit Power Display.....	14
Mechanical and Thermal Indices.....	15
Imaging Functions that Change Acoustic Output	16
Transducer Surface Temperature.....	17
Electrical Safety	18
Level of Protection Against Electrical Shock for Transducers	20
Defibrillators	20
Implantable Devices	20
Possible Combinations with Other Equipment.....	21
Important Information for Maintaining Data Integrity	22
Caring for the Ultrasound System.....	23
Daily Checklist.....	23
Maintenance	24
Repair.....	24
Siemens Healthineers Authorized Care.....	24
Cleaning and Disinfecting the Ultrasound System.....	25
System Surfaces	26
Approved Disinfectant Wipes for the Ultrasound System Surfaces	29
Cleaning the Air Filter.....	32
Caring for Documentation and Storage Devices	33
Caring for Batteries	33
Green Performance Management: Environmental Protection	34
Product Recycling and Disposal	34
Hazardous Substances	34
Recycling Batteries	35
Disposing of the Packaging Materials.....	35
Handling and Disposal of Components and Accessories	36
Energy Efficient Usage: Energy Conservation.....	36

Operating Safety and Environment

Do not operate the ultrasound imaging system until you fully understand the safety considerations and procedures presented in this manual.

Siemens Healthineers recommends the use of gloves during the operation, cleaning, and care of the ultrasound system.

See also: For biohazard considerations related to the use of latex products, refer to page 2-11.

Serious Incident Reporting

If you become aware of any (potential) malfunction of our device which has resulted or could result in serious health consequences for the user, patient or any other person please inform us immediately and if legally required, inform the Competent Authority of your country.


System Symbols

Refer to this table to identify important symbols located on the ultrasound imaging system and transducers.

Symbol	Explanation
	AC (alternating current) voltage source
	Caution: Risk of electric shock (black and yellow illustration)
	Defibrillation-proof type CF applied part
	Type BF applied part
	Consult Instructions for Use
	It is mandatory you refer to the manual. (blue illustration)
	Magnetic resonance (MR) unsafe. Prohibits the use or presence of the medical device in a magnetic resonance imaging (MRI) environment. (black and red illustration)
	Caution (black illustration)
	General warning (black and yellow illustration)
	Attention (black illustration)
	Protective earth ground
	Signal earth ground

Symbol	Explanation
	USB connection
	Ethernet connection
	Continuous wave (pencil) transducer port
	Unlocked position
	Locked position
	Swivel for locking or unlocking the monitor arm
	Partial power on/off control
	ON position of the mains switch (circuit breaker on the ultrasound system)
	OFF position of the mains switch (circuit breaker on the ultrasound system)
	Equipotential terminal Located on the Power Supply Unit (PSU) of the ultrasound system.
	Signal output
	DVI output
	S-Video output
	DVR
STATUS 	Status indicator light for service diagnostics
	Input/output service status
	ECG signal input connection
Aux 1	Auxiliary in/out port
ECG	Electrocardiogram (EKG)
IPX8	Protected against the effects of continuous immersion in water

Symbol	Explanation
	Caution: In the United States of America, federal law restricts this device to sale or use by, or on the order of, a physician.
	Safe working load "XXX" represents the weight of the system unpacked and ready for operation with documentation devices installed Additional objects or devices on the system cannot exceed the safe working load
	Identifies the importer of the medical device
	Authorized representative in Switzerland
	European Conformity mark Indicates the manufacturer declares the product complies with the applicable requirements in Regulation (EU) 2017/745 and other applicable European Union harmonization legislation.
	Manufacturer's declaration of product compliance with applicable EEC directives
	UL symbol for listing as recognized components for Canada and United States of America
	Eurofins MET mark for Canada and United States of America Indicates the product complies with Underwriters Laboratories (UL) and Canadian Standards Association (CSA) safety standards. Includes the following certifications under the mark: E115588 Classified with respect to electrical shock and fire hazards only according to AAMI ES 60601-1 CSA C22.2 No. 60601-1
	UL classified symbol for Canada and United States of America
	Authorized representative in the European community
	EAC (Eurasian Conformity mark) The Eurasian Conformity mark is a product marking to indicate that the product conforms to all Technical Regulations of the Customs Union assessment procedures

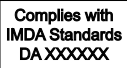
Symbol	Explanation
	Equipment complies with the standard and specifications published by IMDA (Info-communications Media Development Authority of Singapore) "XXXXXX" represents the dealer's license number
	Agência Nacional de Telecomunicações (ANATEL) certification label in Brazil Includes the certification number
	National Communications Commission (NCC) certification of type approval label in Taiwan Includes the approval number
	Technical Conformity Mark for specified radio equipment in Japan
IC	ISED Canada (IC) Innovation, Science and Economic Development (ISED) certification of radio apparatus in Canada Includes the certification number
FCC ID	Federal Communications Commission (FCC) identifier for certification of radio frequency devices in the United States of America Includes the certification number
CMIIT ID	China Ministry of Industry and Information Technology (CMIIT) identifier for radio type certification in China Includes the certification number
IFETEL	Instituto Federal de Telecomunicaciones (IFT) certification in Mexico Includes the certification number
	TÜV Rheinland INMETRO Certification Mark (Brazil)
	TÜV SÜD INMETRO Certification Mark (Brazil)
	UK CA mark Product complies with applicable annexes and directives of the United Kingdom Conformity Assessment
	Ukrainian conformity mark
	The monitor is compliant with the limits for a Class B digital device, pursuant to Part 15 of the Federal Communication Commission (FCC) rules in the United States of America.

Symbol	Explanation
	Product complies with applicable Australian or New Zealand standards for EMC and radio-communications.
	C-Tick Mark Product complies with applicable Australian or New Zealand standards for EMC and radio-communications
	National certification for Korea
IVK	(Installierte Volumen Komponente) Identifier of selected system components or parts for product traceability
System IVK	(System Installierte Volumen Komponente) Identifies the system for product traceability
	Barcode
	Barcode with product identification information according to the following barcode specification: GS1 DataMatrix
	Environmentally friendly use period
	Manufacturer
	Date of manufacture Includes the year and month adjacent to the symbol
	Identifies the manufacturer's serial number of the medical device Includes the number adjacent to the symbol
	Identifies the manufacturer's catalog reference number of the medical device Includes the number adjacent to the symbol
	Indicates the device consists of several components that can be assembled by the manufacturer in multiple configurations
	Indicates the system (a combination of products)
	Indicates the unique device identifier (UDI)
	Indicates the equipment is a medical device
GTIN	Global Trade Item Number
Rev REV	Model revision
Serial SERIAL	Serial number
Count COUNT	Quantity of product
Model MODEL	Product model number


Symbol	Explanation
REF	Device part number
YYYY-MM-DD	Date in YYYY-MM-DD (year-month-day) format
Origin KR ORIGIN KR	Manufacturer country of origin: Republic of Korea
Origin US ORIGIN US	Manufacturer country of origin: United States of America
Origin JP ORIGIN JP	Manufacturer country of origin: Japan
Origin CN ORIGIN CN	Manufacturer country of origin: China
Origin MY ORIGIN MY	Manufacturer country of origin: Malaysia
Origin DO ORIGIN DO	Manufacturer country of origin: Dominican Republic
Origin SK ORIGIN SK	Manufacturer country of origin: Slovakia
MADE IN USA	Location of device manufacturing: United States of America
MADE IN KOREA	Location of device manufacturing: Republic of Korea
MADE IN MALAYSIA	Location of device manufacturing: Malaysia
	Do not install wet
	The wood packaging material for shipping the ultrasound system complies with the International Standards for Phytosanitary Measures No. 15. (example)
	Storage temperature range (example)
	Do not expose to direct sunlight
	Storage atmospheric pressure range (example)
	Storage humidity range (example)
	Indicates this side up
	Do not stack
	Shipping weight (example)

Symbol	Explanation
	Do not allow to get wet
	Fragile. Handle with care.
	Refer to the operator's manual for information about compatible needle guides.
	The product must be properly disposed of in accordance with local, state, and regional laws and regulations. Products bearing this symbol are subject to the European Community Directive 2002/96/EC on waste electrical and electronic equipment (WEEE), amended by Directive 2012/19/EU. For collection and disposal of the product, its components, or its accessories, contact your local customer service representative.
	Control panel height adjustment
	Control panel rotation adjustment
	Do not lean against the cart (black and red illustration)
	Do not exert excessive force when moving or positioning the ultrasound system. (black and red illustration)
	Rear basket weight restriction 5kg MAX 11lb MAX
	Pinch hazard (black and yellow illustration)
	Recycle
	Do not re-use Indicates the medical device is for single use only and must not be used more than one time
	Location of air filter
OFF	Lock released for the wheels
ON	Lock engaged for the wheels
STOP	Lock and swivel lock engaged for the wheels
	Swivel lock engaged for the wheels
	Adjusts the imaging depth or focus
	Adjusts the imaging depth or focus

Labels

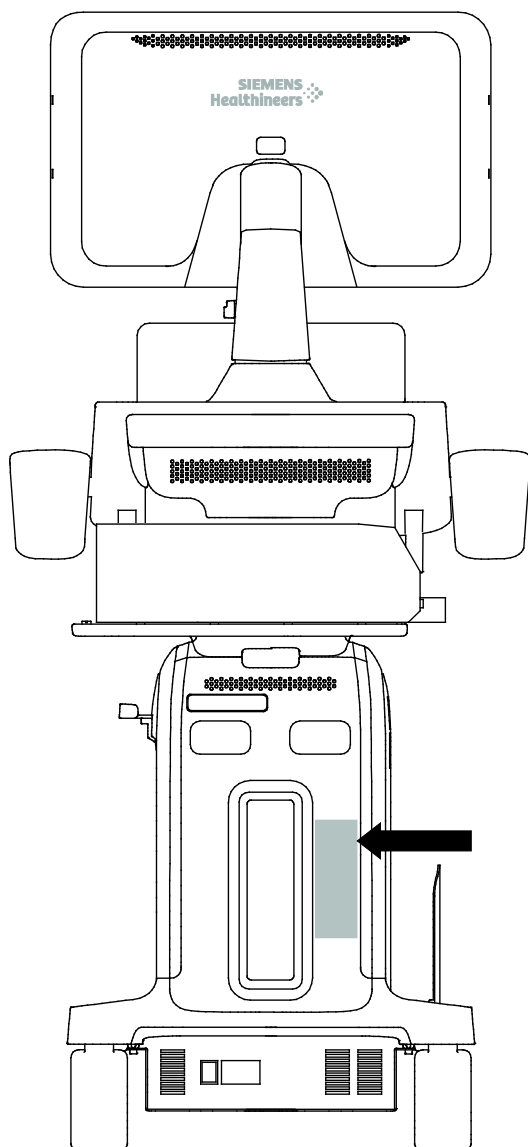


Example of the location of the system label.

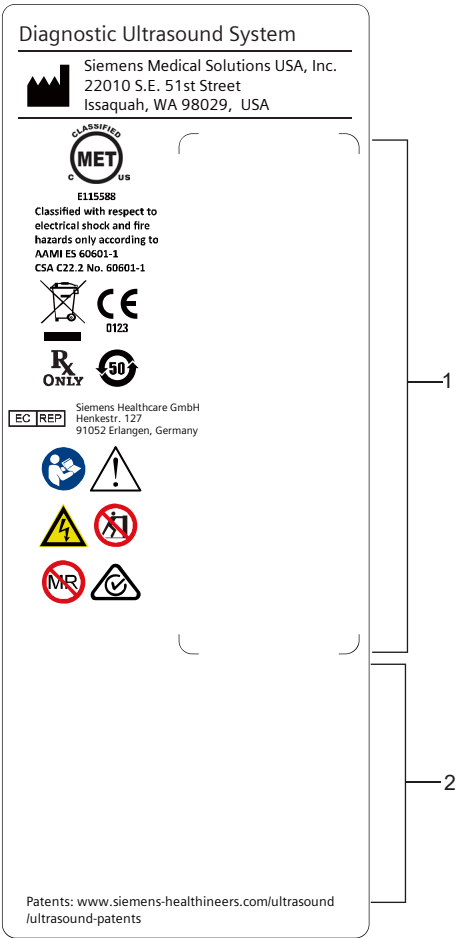


Example of ultrasound system labels with certification labels.

- 1 Identification label
- 2 Country-specific labels

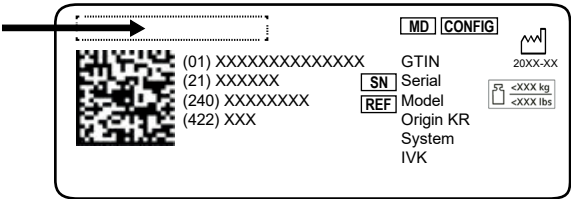

Example of identification label.

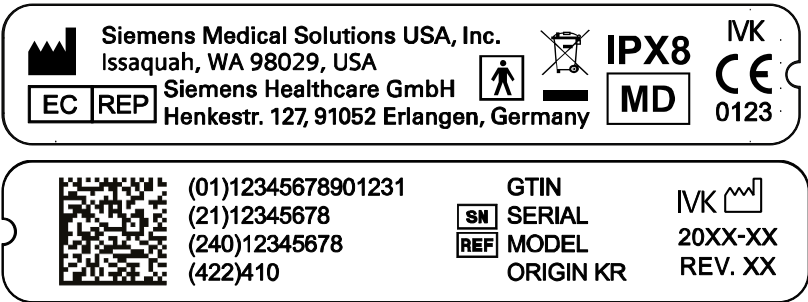

Example of transducer labels.

Biohazard Considerations

-  **WARNING:** This equipment is not suitable for intracardiac use or direct cardiac contact.
-  **WARNING:** For neonatal head imaging, Siemens Healthineers recommends that you exercise special care during neonatal cephalic scanning to avoid possible damage to the posterior region of the eye. The ultrasound energy emitted by the transducer easily penetrates the fontanelles of the infant.
-  **WARNING:** Siemens Healthineers makes every effort to manufacture safe and effective transducers. You must take all necessary precautions to eliminate the possibility of exposing patients, operators, or third parties to hazardous or infectious materials. These precautions should be considered in the use of any application that may indicate the need for such care, and during endocavity or intraoperative scanning; during biopsy or puncture procedures; or when scanning patients with open wounds.
-  **WARNING:** To eliminate the possibility of exposing patients, operators, or third parties to hazardous or infectious materials, always dispose of hazardous or infectious materials according to local, state, and regional laws and regulations.
-  **WARNING:** There have been reports of severe allergic reactions to medical devices containing latex (natural rubber). Health care professionals are advised to identify latex-sensitive patients and to be prepared to treat allergic reactions promptly. For additional information in the U.S.A., refer to FDA Medical Alert MDA91-1.
-  **WARNING:** Ultrasound energy is transmitted more effectively through water than through tissue. When using a standoff device of any kind, for example, a gel pad, the actual mechanical and thermal indices, MI and/or TI, may be higher than indicated in the output display on the system.

The assessment of the biological effects of diagnostic ultrasound on humans is a subject of ongoing scientific research. This system, and all diagnostic ultrasound procedures, should be used for valid reasons, for the shortest possible period of time, and at the lowest mechanical and thermal indices necessary to produce clinically acceptable images.

According to the ALARA (As Low As Reasonably Achievable) principles, acoustic output should be set to the ***lowest level required to satisfactorily perform the examination***.

The ultrasound imaging system complies with the standards of the American Institute of Ultrasound in Medicine (AIUM), the National Electrical Manufacturer's Association (NEMA), the guidelines of the United States Food and Drug Administration (FDA), and the guidelines of the International Electrotechnical Commission (IEC) in terms of safety and acoustic output levels. The ultrasound output levels are stated to permit the user to critically evaluate the system settings in the event of new research findings being announced.

Note on Fetal Examinations

The following recommendation is excerpted from the National Institute of Health in the United States of America. *Consensus Statement on the Use of Ultrasound Imaging During Pregnancy*, Volume 5, No. 1, based on the recommendation issued at the Health Consensus Development Conference, February, 1984:

Ultrasound examination in pregnancy should be performed for a specific medical indication. The data on clinical efficacy and safety do not allow a recommendation for routine scanning at this time.

Ultrasound examination performed solely to satisfy the family's desire to know the fetal sex, to view the fetus, or to obtain a picture of the fetus should be discouraged. In addition, visualization of the fetus solely for educational or commercial demonstrations without medical benefit should not be performed.

In August 1994, the Food and Drug Administration (FDA) notified the medical community and the ultrasound industry regarding its concerns about the misuse of diagnostic ultrasound equipment for non-medical purposes, and to discourage patients from having sonograms for non-medical reasons.


The American Institute of Ultrasound in Medicine (AIUM) has also advocated the responsible use of diagnostic ultrasound for all fetal imaging (August 2005).

Acoustic Output

⚠ WARNING: Ultrasound procedures should be used for valid reasons, for the shortest period of time, and at the lowest mechanical/thermal index setting necessary to produce clinically acceptable images.

Note: For systems distributed in the United States of America, refer to the *Medical Ultrasound Safety* ultrasound education program brochure produced by the AIUM that is shipped with the ultrasound system.

See also: For acoustic output information, refer to Appendix D in this manual.

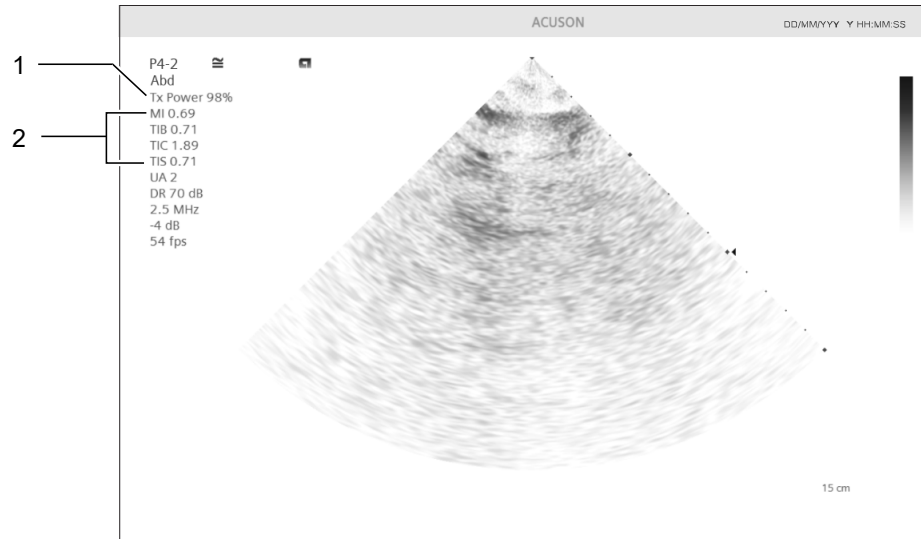

Location of the mechanical and thermal indices and the transmit power display on the image screen.

- 1 Transmit power
- 2 Mechanical and thermal indices

Transmit Power Control

Adjust the transmit power and the corresponding acoustic pressure delivered through the transducer to the patient by using the designated control on the system. It is the main system function that determines the transmitted intensity of ultrasound for all transducers and imaging modes during real-time imaging, though it is not the only function that affects the mechanical and thermal indices. The range and especially the maximum level of the mechanical and thermal indices differ depending on the transducers. In addition, each diagnostic exam type has preset values for mechanical and thermal indices.

See also: Imaging Functions that Change Acoustic Output, page 2-16

Note: Maximum transmit acoustic intensity and the mechanical index for each exam type are limited in accordance with the United States Food and Drug Administration's (FDA) recommendations and guidelines. System default transmit intensity and mechanical index values are always below the FDA recommendations for each exam type. Although some exam types may default to a condition of maximum allowable transmit power, there are other system controls or functions that could raise acoustic output levels.

To increase the transmit power during real-time imaging:

- Rotate **[Tx Power]** clockwise to increase transmit power.

The system briefly highlights the power value to indicate the change.

To decrease the transmit power during real-time imaging:

- Rotate **[Tx Power]** counterclockwise to decrease transmit power.

The system briefly highlights the power value to indicate the change.

Transmit Power Display

The transmit power range is from 0.20% to 100%. Selecting 100%, in combination with other system controls or functions, generates the maximum acoustic intensity and mechanical index for each transducer, where:

$I_{SPTA,3} : \leq 720 \text{ mW/cm}^2$ and $MI \leq 1.9$

Mechanical and Thermal Indices

The ultrasound system incorporates an output display of mechanical and thermal indices to allow you to monitor, and to limit, the amount of ultrasound energy that is transferred to the patient.

The system displays the mechanical and thermal indices during real-time imaging, in all imaging modes, when the mechanical index or the thermal indices are equal to or exceed a value of 0.1.

Use the configuration settings to select the thermal indices for display during neonatal head exams.

Indices display in the abbreviated form shown below:

- **MI:** Mechanical Index
- **TIB:** Bone Thermal Index (fetal application)
- **TIS:** Soft Tissue Thermal Index
- **TIC:** Cranial Bone Thermal Index

Imaging Functions that Change Acoustic Output

 **WARNING:** Observe the real-time display of mechanical and thermal indices (MI/TI) at all times.

In addition to the adjustment of the transmit power, adjustment of the following imaging functions and/or controls may affect the acoustic output:

- Automatic Time-out
- M-mode cursor
- Color and Power ROI Position and Size; Steering Angle for linear array transducers
- Doppler Gate Position and Size; Steering Angle for linear array transducers
- Doppler PRF
- Color PRF
- Exam Type
- Field of View (Scan Angle); 2D Steering Angle for linear array transducers
- Focus
- Zoom
- Frame Rate
- Freeze
- Image Depth
- Imaging Mode
- Multi-Frequency
- Power On/Off
- Image presets
- Harmonics
- Line Density
- Simultaneous/Update
- Transducer

Transducer Surface Temperature

The following table provides the maximum surface temperature of the transducers compatible with the ultrasound system.

The maximum surface temperature of the transducers is in accordance with IEC 60601-2-37.

Transducer	Maximum Temperature	
	TMM	Still Air
5C1	42.6°C	37.4°C
6C1	41.3°C	39.2°C
7C2	42.5°C	40.9°C
12L3	38.9°C	42.1°C
11L4	39.7°C	38.7°C
13L4	41.7°C	36.7°C
14L4	42.7°C	41.4°C
16L4	40.0°C	39.4°C
18H5	39.8°C	41.4°C
11M3	41.1°C	40.5°C
L10-5v	39.7°C	37.1°C
5P1	42.0°C	44.7°C
8V4	41.3°C	39.5°C
10V4	41.7°C	40.0°C
5VT	41.5°C	31.5°C
P4-2	42.5°C	39.0°C
9MC3	42.3°C	40.3°C
10MC3	42.5°C	39.7°C
10BCC3	40.9°C	34.4°C
9VC2	38.3°C	37.8°C
8VC3	40.1°C	43.5°C
9VE3	42.5°C	37.2°C
9VE4	40.1°C	36.8°C
CW2	35.3°C	40.1°C
CW5	34.8°C	31.9°C
CW8	38.3°C	31.2°C

TMM = Tissue Mimicking Material

Electrical Safety

-  **WARNING:** To avoid electric shock, use a protective earth connection to connect the ultrasound system to the power supply (supply mains). The protective earth connection ensures that the circuit breaker on the ultrasound system (mains switch) will disconnect the power supply in the event of a short circuit.
-  **WARNING:** For systems used in the U.S.A.: To ensure grounding reliability, only plug the system into a hospital-grade power outlet (mains outlet).
-  **WARNING:** The power connector (mains plug) for the ultrasound system is a three-prong grounded plug in the U.S.A. and should never be adapted to any two-prong, non-grounded outlet, either by modifying the plug or by using an adapter. In the U.S.A., proper grounding requires the power connector to be plugged into a hospital-grade power outlet (mains outlet).
-  **WARNING:** To avoid electrical shock, never modify the ultrasound system's power connector (mains plug), as doing so may overload your facility's power circuits. To ensure grounding reliability, connect the system only to an equivalent outlet.
-  **WARNING:** To avoid electrical shock, never use equipment or a mains power cord that shows signs of wear or tampering, or whose ground plug has been bypassed using an adapter.
-  **WARNING:** Equipment connected to the ultrasound system and in the patient zone must be powered from a medically-isolated power source or must be a medically-isolated device. Equipment powered from a non-isolated source can result in chassis leakage currents exceeding safe levels. Chassis leakage current created by an accessory or device connected to a non-isolated outlet may add to the chassis leakage current of the ultrasound system.
-  **WARNING:** Using an extension cord or multi-socket outlet setup to provide power to the ultrasound system, or to the system's peripheral devices may compromise the system grounding and cause your system to exceed leakage current limits.
-  **WARNING:** To avoid electrical shock and damage to the ultrasound system, power off and unplug the equipment from the power supply (supply mains) before cleaning and disinfecting the system.
-  **WARNING:** To avoid electrical shock and damage to the control panel resulting from ingress of liquid, place the gel and gel warmer on the side of the system closest to the patient.
-  **WARNING:** Do not pour any fluid onto the ultrasound system surfaces, as fluid seepage into the electrical circuitry may cause excessive leakage current or system failure.

-  **WARNING:** To ensure proper grounding and leakage current levels, it is the policy of Siemens Healthineers to have an authorized Siemens Healthineers representative or Siemens Healthineers-approved third party perform all on-board connections of documentation and storage devices to the ultrasound system.
-  **WARNING:** To maintain the safety and functionality of the ultrasound system, electrical safety tests must be performed as specified by local safety regulations, or as needed.
-  **WARNING:** To prevent excessive leakage current from contacting the patient, do not touch a user-accessible connector on the system while touching or scanning the patient. User-accessible connectors include the ECG connector, a transducer connector, a USB connector, and any other audio, video, or data transmission connectors.
-  **WARNING:** Connecting peripheral devices to accessory outlets on the ultrasound system effectively creates a medical electrical system, resulting in a reduced level of safety.
-  **WARNING:** Do not modify this equipment without authorization from Siemens Healthineers.
-  **Caution:** To avoid the possibility of static shock and damage to the system, avoid the use of aerosol spray cleaners on the monitor screens.
-  **Caution:** Do not use spray cleaners on the ultrasound system, as this may force cleaning fluid into the system and damage electronic components. It is also possible for the solvent fumes to build up and form flammable gases or damage internal components.
-  **Caution:** To maintain the safety and functionality of the ultrasound system, maintenance must be performed every 24 months. Electrical safety tests must also be performed at regular intervals as specified by local safety regulations, or as needed.

Level of Protection Against Electrical Shock for Transducers

⚠ WARNING: Only use Type BF transducers with the ultrasound system to maintain a level of protection against electrical shock.

According to IEC 60601-1, the assemblies for the linear, curved, and phased array transducers provide a "Level of Protection Against Electrical Shock" of "Type BF."



The Type BF icon is located on the transducer label.

Defibrillators

⚠ WARNING: The ECG function is designed to withstand the effects of defibrillation. However, when possible, disconnect the ECG leads during defibrillation since a malfunction of the safety controls could otherwise result in electrical burns for the patient.

For patient safety, be sure to use defibrillators that do not have grounded patient circuits.



The Type CF defibrillation proof icon is located on the ECG module.

Implantable Devices

⚠ WARNING: Ultrasound systems, like other medical equipment, use high-frequency electrical signals that can interfere with implantable devices such as pacemakers and implantable cardioverter-defibrillators (ICDs). If the patient has such an implantable device, you should be aware of any interference in its operation and immediately power off the ultrasound system.

Possible Combinations with Other Equipment

⚠ WARNING: Accessory equipment connected to the analog and digital interfaces must be certified according to the respective EN or IEC standards (for example, EN 60950 or IEC 60950 for data processing equipment and EN 60601-1 or IEC 60601-1 for medical equipment). Anyone who connects additional equipment to any of the signal input or signal output ports configures a medical system and is therefore responsible that the system complies with the requirements of the system standards EN 60601-1 or IEC 60601-1. Siemens Healthineers can only guarantee the performance and safety of the devices listed in the *Instructions for Use*. If in doubt, consult the Siemens Healthineers service department or your local Siemens Healthineers representative.

The ultrasound system supports three on-board documentation devices. The ultrasound system can support two on-board documentation devices at one time.

On-board peripheral devices must be installed by an authorized Siemens Healthineers representative or by a Siemens Healthineers-approved third party. Devices installed by other people will be at the user's risk and may void the system warranty.

Barcode Reader

⚠ WARNING: To avoid electrical shock, do not connect the barcode reader to the system when more than two on-board documentation devices are connected to the system.

⚠ WARNING: Do not stare into the beam of the barcode reader. This action can create a risk for injury to the user.

You can connect a Siemens Healthineers-approved barcode reader directly to a USB port on the ultrasound system. The barcode reader provides direct entry of patient information in the registration form, for example, patient name, patient identifier, or physician identifier.

Important Information for Maintaining Data Integrity

 **WARNING:** To prevent the loss of data that results from unplugging the system from the supply mains (power source) before the system has completed its shutdown sequence, follow the system shutdown instructions in this manual.

 **Caution:** To prevent corruption of data stored on inserted media, do not interrupt the exporting process. Do not delete exported patient data from the local database until the exporting process is complete and the inserted media is determined to be full (or no more exporting is planned for the inserted media).

- To prevent the loss of data that results from power failures and other system "down" occurrences, you must archive important data, such as patient records, onto an external media, for example, a USB storage device.
- Loss of data is to be expected and its retrieval is not normally possible under the following conditions: loss of power to the ultrasound system, hard disk failure, CPU failure, system lockup, and other similar causes.
- Should an abnormal system shutdown occur, retrieval of data not saved to the hard disk or not archived to an external recording media is not normally possible.

An abnormal system shutdown occurs if you do not power off the ultrasound system using the partial power on/off control located on the control panel. Other examples of abnormal system shutdown include: equipment malfunction, loss of power, or pressing and holding the partial power on/off control longer than four seconds.

- Should an abnormal system shutdown occur, the system may initially require additional time to reboot or to respond to user input. This is due to the operating system performing a background scan of the hard disk to detect and segregate any truncated or corrupted files.

Caring for the Ultrasound System

It is the responsibility of the user to verify that the ultrasound system is safe for diagnostic operation on a daily basis. Each day, prior to using the system, perform each of the steps in the Daily Checklist.

All exterior parts of the system, including the control panel, keyboard, transducers, and biopsy devices should be cleaned and/or disinfected as necessary or between uses. Clean each component to remove any surface particles. Disinfect components to destroy pathogens and other microorganisms.

See also: Cleaning the Air Filter, page 2-32

Daily Checklist

⚠ WARNING: To minimize the risk of cross-contamination and infectious diseases, a sterile, non-pyrogenic transducer sheath must be in place during procedures requiring sterility.

⚠ WARNING: To avoid electrical shock, you must visually inspect a transducer prior to use. Do not use a transducer with cracked, punctured, or discolored casing or a frayed cable.

Discoloration Exception: The use of approved cleaners and disinfectants may cause discoloration of the transducer. You can continue to use a transducer discolored due to the use of approved cleaners and disinfectants only.

See also: For additional information about a list of disinfectants for transducer, refer to Chapter 6 in this manual.

Perform the following each day before using the ultrasound system:

- Visually inspect all transducers. Do not use a transducer with a cracked, punctured, or discolored casing or a frayed cable. You can continue to use a transducer discolored due to the use of approved cleaners and disinfectants only.

Verify the electrical component of the transducer connector is clean and free from gel or particles.

- Visually inspect all power cords. Do not turn on the power if a cord is frayed or split, or shows signs of wear.

If your system's power cord is frayed or split, or shows signs of wear, contact your Siemens Healthineers service representative for power cord replacement.

- Visually inspect the ECG connector and the cable. Do not use the ECG function if the connector or cable is damaged or broken.
- Verify the trackball, DGC slide controls, and other controls on the control panel are clean and free from gel or other contaminants.

Once the system is powered on:

- Visually check the on-screen displays and lighting.
- Verify the screen displays the current date and time.
- Verify the transducer identification and indicated frequency are correct for the active transducer.

Maintenance

 **Caution:** To maintain the safety and functionality of the ultrasound system electrical safety tests must be performed as specified by local safety regulations, or as needed.

Repair

 **WARNING:** Do not modify this equipment without authorization from Siemens Healthineers.

For questions regarding repair or replacement of any equipment parts on your system, contact your Siemens Healthineers service representative.

Siemens Healthineers Authorized Care

Installers and operators must observe any statutory regulations that govern the installation, operation, inspection, and maintenance of this equipment.

To ensure the safety of patients, operators, and third parties, the equipment must be inspected, every 24 months and the replacement of parts is performed as necessary. This maintenance must be performed by a qualified Siemens Healthineers authorized representative. It is important to inspect the equipment more frequently if it is operated under extraordinary conditions.

Perform inspections and maintenance at the prescribed intervals to avoid worn and hazardous parts due to wear. Contact the Siemens Healthineers service department for the required maintenance.

As manufacturers and installers of ultrasound equipment, Siemens Healthineers cannot assume responsibility for the safety properties, reliability, and/or performance of the equipment, if:

- Installations, extensions, readjustments, modifications, additions, or repairs are carried out by persons not specifically authorized by Siemens Healthineers.
- Components that affect the safe operation of the system are replaced by parts not authorized by Siemens Healthineers.
- The electrical installation of the room where the equipment is located does not meet the power and environmental requirements stated in this manual.
- The equipment is not used in accordance with the operating instructions.
- The system is operated by personnel not adequately educated or trained.

Siemens Healthineers suggests that you request any person who performs maintenance, or repairs, to provide you with a certificate showing:

- The nature and extent of the work performed
- Changes in rated performance
- Changes in working ranges
- Date of service
- Name of person or firm performing the service
- Signature of person performing the service

Cleaning and Disinfecting the Ultrasound System

You must take all necessary precautions to eliminate the possibility of exposing patients, operators, or third parties to hazardous or infectious materials. Use universal precautions when cleaning and disinfecting. You should treat all portions of the imaging system that come in contact with human blood or other body fluids as if they were known to be infectious.

All exterior parts of the system, including the control panel, keyboard, transducers, and biopsy devices should be cleaned and/or disinfected as necessary or between uses. Clean each component to remove any surface particles. Disinfect the components to destroy pathogens and other microorganisms.

-  **WARNING:** To avoid electric shock and damage to the system, always power off and disconnect the equipment from the AC power source before cleaning and disinfecting.
-  **WARNING:** Disinfectants and cleaning methods listed are recommended by Siemens Healthineers for compatibility with product materials, not for biological effectiveness. Refer to disinfectant label instructions for guidance on disinfection efficacy and appropriate clinical uses.
-  **WARNING:** The use of any disinfectants other than those specified in the instructions for use may damage the ultrasound system and accessory surfaces and, as a result, may create electrical hazards for the patients and/or users.
-  **WARNING:** Contents of some disinfecting agents are known to be health hazards. Their concentration in the air must not exceed an applicable specified limit. Comply with the manufacturer's instructions when using these agents.
-  **Caution:** Do not clean or disinfect the system with chlorinated or aromatic solvents, acidic or basic solutions, isopropyl alcohol or strong cleaners such as ammoniated products, as these can damage the surface of the system. Use the recommended cleaning and disinfecting procedure. Only use approved cleaners and disinfectants on the ultrasound system.
-  **Caution:** Do not use spray cleaners on the ultrasound system as this may force cleaning fluid into the system and damage electronic components. It is also possible for the solvent fumes to build up and form flammable gases or damage internal components.
-  **Caution:** Do not pour any fluid onto the ultrasound system surfaces, including the control panel and touch screen, as fluid seepage into the electrical circuitry may cause excessive leakage current or system failure.
-  **Caution:** Prior to cleaning or disinfecting the system, you must disconnect and remove transducers from the ultrasound system. Unintended contact with disinfectants approved only for use with the system can result in damage to the transducers.



Note: Iodine and betadine may discolor the control panel. There is no associated degradation of imaging performance or control panel reliability due to discoloration by iodine or betadine.

System Surfaces

See also: For instructions on removing and installing transducer holders and the gel warmer, refer to Chapter 3 in this manual.

To	Do This
Clean and disinfect the surfaces of the ultrasound system	 Caution: To avoid the possibility of static shock and damage to the ultrasound system, avoid the use of aerosol spray cleaners on the monitor screen and touch screen. 1. Power off the ultrasound system and unplug the power cord from the power supply (supply mains). 2. Use a clean gauze pad or lint-free cloth, lightly moistened with a mild detergent, to wipe the surfaces of the ultrasound system. – Continue wiping the surfaces of the ultrasound system until there is no visible residue. – Gently wipe the surface of the touch screen and monitor. – Take particular care to clean the areas near the trackball and the slide controls. Ensure these areas are free of gel and any other visible residue. – Ensure that fluids do not seep into any openings on the ultrasound system, such as the seams on the control panel. 3. After cleaning, use a clean, lint-free cloth to dry the surfaces. 4. Use an approved disinfectant wipe to disinfect the ultrasound system and accessories, except the transducer ports. 5. After cleaning and disinfecting, reconnect the ultrasound system power cord to the power supply (supply mains).
Clean the array transducer ports on the ultrasound system	 Caution: To avoid damage to the transducer ports, never use gauze pads, lint-free cloths, or solutions to clean the transducer ports. 1. Disconnect transducers from the array transducer ports. 2. Use compressed air to remove particles from the ports. 3. Reconnect the transducer and then verify the transducer identification and indicated frequency are correct for the active transducer.
Clean and disinfect the holders for transducers and gel	1. Remove the holder from the ultrasound system. 2. Clean the holder under running water, using a mild detergent and then dry with a clean, lint-free cloth. 3. Use an approved disinfectant wipe to disinfect the holder. 4. Reattach the holder to the ultrasound system.

To	Do This
Clean and disinfect the rear basket or storage bin	See also: For information about removing and attaching the rear basket or storage bin, refer to Chapter 3 in this manual 1. Remove the rear basket or storage bin from the ultrasound system. 2. Use a clean gauze pad or lint-free cloth, lightly moistened with a mild detergent, to wipe the basket and the point of attachment on the ultrasound system. 3. After cleaning, use a clean, lint-free cloth to dry the basket or storage bin and the point of attachment on the ultrasound system. 4. As required, use an approved disinfectant wipe to disinfect the basket and the point of attachment on the ultrasound system. 5. Reattach the rear basket or storage bin to the ultrasound system.
Clean and disinfect the gel warmer	 Caution: To avoid damage to the gel warmer, do not immerse the gel warmer in water or any solution. Note: To facilitate cleaning, remove the cap at the bottom of the gel warmer. Rotate the cap counterclockwise to remove the cap; rotate the cap clockwise to secure the cap. 1. Remove the gel warmer from the ultrasound system. 2. Unplug the power cable from the gel warmer. 3. Clean and disinfect the gel warmer using an approved disinfectant wipe. 4. Reattach the gel warmer to the ultrasound system.

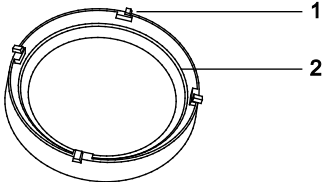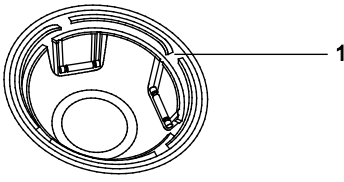To	Do This
Clean and disinfect the trackball	⚠ Caution: Do not drop or place foreign objects inside the trackball assembly because doing so may affect the trackball operation and damage the ultrasound system.  <i>Example of the trackball ring.</i> 1 Tab for inserting into the trackball assembly 2 Raised markings on either side of the ring  <i>Example trackball assembly.</i> 1 Point of attachment for the trackball ring 1. Rotate the ring around the trackball counterclockwise and then carefully lift up the ring to remove the ring and trackball. 2. Clean the ring and trackball with a cotton swab or lint-free pad moistened with mild detergent solution. 3. Clean the inside of the trackball assembly with a cotton swab or lint-free pad moistened with mild detergent solution. 4. As required, use an approved disinfectant wipe to disinfect the ring, trackball, and trackball assembly. 5. Allow the trackball components to completely dry before reassembly. 6. Reinstall the ring and trackball. a. Place the trackball inside the trackball assembly. b. Place the ring over the trackball, taking care to align the tab with the point of attachment on the trackball assembly. c. Rotate the ring clockwise until the ring snaps into place and the ridges of the ring are in a horizontal position.

Approved Disinfectant Wipes for the Ultrasound System Surfaces

Use only approved disinfectant wipes on the surfaces of the ultrasound system and components, except for the monitor screen, touch screen, and transducer ports.

Note: For spray disinfectants, you must apply the disinfectant to a clean, lint-free cloth before wiping the surfaces of the ultrasound system.

	Monitor Screen	Touch Screen	Trackball Assembly	Gel Warmer	All other ultrasound system surfaces*
Accel INTERVention Wipes	✓	✓	✓	✓	✓
Acrylan	✓	✓	✓	✓	✓
ANIOS OXI'ZYME	✓	✓	✓	✓	✓
ANIOXY-TWIN	✓	✓	✓	✓	✓
ANIOXYDE 1000	✓	✓	✓	✓	✓
ANIOXYDE 1000 LD	✓	✓	✓	✓	✓
antifect N liquid	✓	✓	✓	✓	✓
AZO 70% IPA Wipes & Spray	✓	✓	✓	✓	✓
AzoMax wipes & spray	✓	✓	✓	✓	✓
Azo Universal Cleaning and Disinfectant Wipes	✓	✓	✓	✓	✓
Bacillol 30 Tissues & Foam	✓	✓	✓	✓	✓
Bacillol AF	✓	✓	✓	✓	✓
Bacillol Tissues	✓	✓	✓	✓	✓
Bacillol Wipes	✓	✓	✓	✓	✓
CaviCide	✓	✓	✓	✓	✓
CaviCide1	✓	✓	✓	✓	✓
CaviWipes	✓	✓	✓	✓	✓
CaviWipes 1	✓	✓	✓	✓	✓
CaviWipes 2.0	✓	✓	✓	✓	✓
CLEANISEPT WIPES	✓	✓	✓	✓	✓
CLEANISEPT WIPES forte	✓	✓	✓	✓	✓
Clinell Clorox Wipes	✓	✓	✓	✓	✓
Clinell Universal Wipes & Spray	✓	✓	✓	✓	✓
Clorox Healthcare Bleach Germicidal Cleaner	✓	✓	✓	✓	✓
Clorox Healthcare Bleach Germicidal Wipes	✓	✓	✓	✓	✓
Clorox Healthcare Fuzion Cleaner Disinfectant	✓	✓	✓	✓	✓
Clorox Healthcare Hydrogen Peroxide Cleaner Disinfectant	✓	✓	✓	✓	✓
Clorox Healthcare VersaSure Cleaner Disinfectant Wipes	NA	✓	✓	NA	✓

	Monitor Screen	Touch Screen	Trackball Assembly	Gel Warmer	All other ultrasound system surfaces*
Dispatch	✓	✓	✓	✓	✓
Ecolab Quaternary Disinfectant Wipes	✓	✓	✓	✓	✓
ESCULASE	✓	✓	✓	✓	✓
gigasept PAA concentrate	✓	✓	✓	✓	✓
gigasept pearls	✓	✓	✓	✓	✓
gigasept powerTRIO	✓	✓	✓	✓	✓
HospiCare 3525S Disinfectant Wipes	✓	✓	✓	✓	✓
HospiCare 60R Multi-surface Wipes	✓	✓	✓	✓	✓
Incidin OxyFoam S	✓	✓	✓	✓	✓
Incidin OxyWipe	✓	✓	✓	✓	✓
Incidin OxyWipe S	✓	✓	✓	✓	✓
INTERCEPT Detergent	✓	✓	✓	✓	✓
Medi-Fect Wipes	✓	✓	✓	✓	✓
Meliseptol Wipes sensitive	✓	✓	✓	✓	✓
Micro-Kill Bleach Germicidal Bleach Wipes	✓	✓	✓	✓	✓
Mikrobac Tissues	✓	✓	✓	✓	✓
mikrozyd AF wipes & liquid	✓	✓	✓	✓	✓
mikrozyd PAA wipes	✓	✓	✓	✓	✓
mikrozyd sensitive wipes & liquid	✓	✓	✓	✓	✓
mikrozyd universal wipes & liquid	✓	✓	✓	✓	✓
mikrozyd universal wipes green line	✓	✓	✓	✓	✓
mycolio lite Alcohol Wipes	✓	✓	✓	✓	✓
Oxivir Tb	✓	✓	✓	✓	✓
PDI Easy Screen Cleaning Wipe	✓	✓	✓	✓	✓
PHAGO'SPORE & PHAGO'WIPES SPORE	✓	✓	✓	✓	✓
PROSAT Sterile Meltblown Polypropylene Wipes	✓	✓	✓	✓	✓
Protex Disinfectant Spray and Wipes	✓	✓	✓	✓	✓
Protex ULTRA Disinfectant Wipes	✓	✓	✓	✓	✓
RAPICIDE PA	✓	✓	✓	✓	✓
Rely+On PeraSafe	✓	✓	✓	✓	✓
Reynard IPA Surface Disinfection Wipes	✓	✓	✓	✓	✓
Reynard Premier Detergent & Disinfectant Wipes	✓	✓	✓	✓	✓

	Monitor Screen	Touch Screen	Trackball Assembly	Gel Warmer	All other ultrasound system surfaces*
Sani-Cloth 70% Wipes	✓	✓	✓	✓	✓
Sani-Cloth Active Wipes	✓	✓	✓	✓	✓
Sani-Cloth AF Universal Wipes	✓	✓	✓	✓	✓
Sani-Cloth AF3 Germicidal Disposable Wipes	✓	✓	✓	✓	✓
Sani-Cloth Bleach Germicidal Disposable Wipe	✓	✓	✓	✓	✓
Sani-Cloth Chlor Wipes	✓	✓	✓	✓	✓
Sani-Cloth Plus Germicidal Disposable Cloth	✓	✓	✓	✓	✓
Sani-Cloth Prime Germicidal Disposable Wipes	✓	✓	✓	✓	✓
Sani Professional No-Rinse Sanitizing Multi-Surface Wipes	✓	✓	✓	✓	✓
Septiwipes	✓	✓	✓	✓	✓
SideKick Disinfecting Wipes	✓	✓	✓	✓	✓
SONO Ultrasound Wipes	✓	✓	✓	✓	✓
STERIMED 100	✓	✓	✓	✓	✓
STERICID Wipes	✓	✓	✓	✓	✓
Super Sani-Cloth Germicidal Disposable Wipe	✓	✓	✓	✓	✓
SURFA'SAFE Premium	✓	✓	✓	✓	✓
terralin PAA	✓	✓	✓	✓	✓
Transeptic	✓	✓	✓	✓	✓
Tristel DUO	✓	✓	✓	✓	✓
Tristel Duo ULT	✓	✓	✓	✓	✓
Tristel Solo for Ultrasound	✓	✓	✓	✓	✓
Tristel ULT	✓	✓	✓	✓	✓
trophon Companion Cleaning Wipes	✓	✓	✓	✓	✓
Tuffie 5 Universal Sanitising Wipes	NA	✓	NA	NA	NA
Virex II 256	✓	✓	✓	✓	✓
Virusolve+ Disinfectant Wipes	✓	✓	✓	✓	✓
Wet Wipe Ethanol Disinfection with detergent	✓	✓	✓	✓	✓
WIP'ANIOS Excel	✓	✓	✓	✓	✓
WipeClean Ethanol Disinfection 80%	✓	✓	✓	✓	✓


*All other ultrasound system surfaces include the surfaces not listed in the above table, including the storage bin and transducer holders.

✓ = approved

NA = not approved

Cleaning the Air Filter

The ultrasound system has a removable, washable air filter. Check and clean the air filters weekly, or as needed, to maintain proper system cooling.

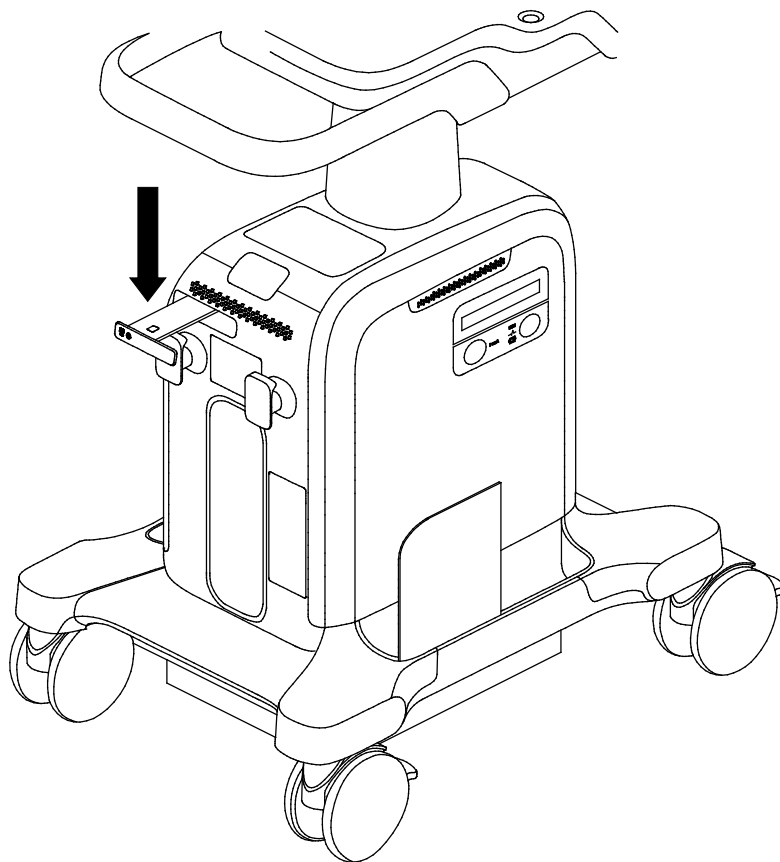

Example of the location of the air filter.

To remove and clean the air filter:

-  **Caution:** Do not scrub, stretch, or bend the filter, or apply heat to the filter, as doing so could damage the filter.
1. Power off the ultrasound system and unplug the power cord from the power supply (supply mains).
 2. Push the air filter tray in to release the locking mechanism and then pull the tray from the system.
 3. Rinse the air filter tray with running water and allow the filter to completely dry.
To hasten drying, you may gently shake the air filter tray, or blot the filter with a clean, lint-free cloth.
-  **Caution:** Do not insert the air filter tray into the ultrasound system with a wet filter as this can damage the ultrasound system.
4. Slide the dry air filter tray back into the ultrasound system.
 5. Plug the power cord into the power supply (supply mains).

Caring for Documentation and Storage Devices

-  **Caution:** Computer viruses on a USB-compatible device can infect the ultrasound system. Check for viruses before connecting the device to the system.

Note: Study data stored directly to a USB device can be lost. Do not use a USB device for permanent data storage.

For information on the care of an optional documentation or storage device, refer to the manufacturer's operating instructions that accompanied the device.

Caring for Batteries

-  **WARNING:** Do not strike or drop batteries, allow batteries to contact water or other fluids, disassemble batteries, allow conductive object to contact a battery's terminals, cause a battery to become short circuited, heat batteries, or expose batteries to fire. Any of these actions can compromise the structural integrity of a battery. Compromising the structural integrity of a battery can result in battery leakage, heat generation, fire, or explosion, causing possible personal injury.
-  **WARNING:** Do not use a battery if it leaks fluid or has changed shape. If skin or clothing comes in contact with fluid from the battery, thoroughly wash the area immediately with clean water. If any fluid comes in contact with a user's eyes, immediately flush their eyes with water and seek medical attention.
-  **WARNING:** Replace batteries with the same or equivalent type. Use of incompatible batteries can result in battery leakage, heat generation, fire, or explosion, causing possible personal injury.

Green Performance Management: Environmental Protection

Siemens Healthineers recognizes its responsibility to minimize the environmental impacts of services, operations, and products. For more information, contact your local Siemens Healthineers representative.

Product Recycling and Disposal

 **WARNING:** Do not attempt to recharge non-rechargeable batteries, such as the batteries included on printed circuit boards. Recycle non-rechargeable batteries according to local, state, and regional laws and regulations.

Dispose of this product according to local, state, and regional laws and regulations.

Batteries and electrical and electronic equipment can contain hazardous substances. If released, the hazardous substances can harm people and the environment.

Siemens Healthineers provides disassembly instructions to treatment facilities for the safe and proper removal and recycling of electrical and electronic components in this product. For more information, contact your local Siemens Healthineers representative.

To the extent required by local, state, and regional laws and regulations, Siemens Healthineers has programs for the return of used products. For more information, contact your local Siemens Healthineers representative.

Hazardous Substances

The 54.5 cm (21.5 inch) flat panel display monitor does not contain mercury.

Recycling Batteries

 **WARNING:** Never dispose of batteries by burning or by flushing into any waste water system, for example, a lavatory. Compromising the structural integrity of a battery can result in leakage or explosion and the potential for personal injury.

 **WARNING:** Do not throw batteries into the trash. Collect and recycle used batteries separate from other waste.

Recycle batteries according to local, state, and regional laws and regulations. Use a battery collection program available in your country to recycle batteries.

To the extent required by local, state, and regional laws and regulations, Siemens Healthineers will collect and recycle batteries for this product at no charge. Contact your local Siemens Healthineers representative for battery shipment instructions.

Disposing of the Packaging Materials

Dispose of or recycle the packaging materials according to local, state, and regional laws and regulations.

To the extent required by local, state, and regional laws and regulations, Siemens Healthineers will collect and dispose of packaging materials for this product. For more information, contact your local Siemens Healthineers representative.

Handling and Disposal of Components and Accessories

⚠ WARNING: Observe local, state, and regional laws and regulations for the disposal of the ultrasound system components and accessories.

⚠ WARNING: Bodily fluids on used needles and needle guides can transmit infectious diseases. To eliminate the possibility of exposing patients, operators, or third parties to hazardous or infectious materials, always dispose of the needle and the needle guide according to local, state, and regional laws and regulations.

Component or Accessory	Handling Instructions
Bracket for needle guides	(Available only for the reusable bracket) Replace the bracket if the bracket is damaged in any way. Refer to the in-box instructions for disposal guidelines for the reusable bracket.
Needle guides	The sterile-packaged needle guide accessories are single-use items. Refer to the in-box instructions for disposal guidelines for the disposable needle guide accessories.
Protective overlay for the control panel	The protective overlay for the control panel is disposable. Replace the protective overlay if the protective overlay is damaged in any way. For example, if the protective overlay has a hole or does not fit on the control panel, replace the overlay.
Air filter tray	See also: Cleaning the Air Filter, page 2-32 The air filter tray is reusable. Replace the air filter tray if the air filter is damaged in any way. For example, the air filter has a hole, or the air filter tray will not fit in the air filter tray slot.

Energy Efficient Usage: Energy Conservation

See also: For information about the QuikStart feature and supplying power to the ultrasound system, refer to Chapter 3 in this manual.

For moderate energy conservation when the ultrasound system is not in use, activate system standby.

For improved energy conservation when the ultrasound system is not in use, power off the system. Keep the system plugged into the power supply (supply mains). Ensure the circuit breaker (mains switch) on the ultrasound system is in the on position (I = ON).

For maximum energy conservation when the ultrasound system is in storage, power off and unplug the system from the power supply (supply mains). Ensure the circuit breaker (mains switch) on the ultrasound system is in the off position (O = OFF). Implementing maximum energy conservation for long periods of time can shorten the life of the system batteries.

3 System Setup

Initial Setup	3
System Ergonomics	3
System Startup	4
Plugging in the Ultrasound System.....	4
Supplying Power to the Ultrasound System.....	5
Activating System Standby	6
Activating the Battery Support Feature.....	8
Recalibrating the Battery Pack	11
Connecting and Disconnecting Transducers	12
ACUSON Juniper System	12
ACUSON Juniper Select System	13
Array Transducers.....	14
Continuous Wave Transducers	15
Protective Transducer Holders and Gel Holders	16
Transducer Cable Management	17
Gel Warmer	18
Installing the Storage Container	20
Adjusting the Audio Volume.....	21
Adjusting the Screen Brightness.....	22
Input/Output Connections	24
Connecting Peripheral Equipment and Accessories	25
Installing the Protective Overlay	28
Connecting the Ultrasound System to a Network	29
Wired Connection.....	29
Wireless Network Connections.....	29
Use of Wireless Adapters	30
Connecting to Remote Services.....	30

Moving the System	31
Securing the Wheel Locks.....	32
Wheels without Swivel Locks	32
Wheels with Swivel Locks	33
Before the Move	34
ACUSON Juniper System	35
ACUSON Juniper Select System	36
During the Move	37
After the Move	37
Shipping the Ultrasound System	38
Installing Updates	39
Installing System Software	40
Installing a Software-Based Option	41
Accessing Help	42

Initial Setup

 **WARNING:** Do not tip the system in any direction more than 10 degrees. This action can cause the system to tip over and create a risk for injury to the user or patient and damage to the system.

 **WARNING:** Do not lean on or apply excessive force to the control panel or monitor. These actions can cause the ultrasound system to tip over and create a risk for injury to the user or patient and damage to the system.

The ultrasound system is initially unpacked and installed by a Siemens Healthineers representative who will verify the operation of the system. Any transducers, documentation and storage devices, accessories, and options delivered with your system are also connected and installed for you.

System Ergonomics

The relationship of repetitive ultrasound scanning with musculoskeletal disorders is indicated in literature for users conducting all types of ultrasound studies. Examples of musculoskeletal disorders include pain and injury to the neck, shoulder, back, wrist, hands, and fingers.

We recommend keeping informed about the trends in literature for continued education regarding scanning techniques, for example, handling a transducer to obtain a quality image and the proper posture to use while scanning a patient.

The ultrasound system provides features designed for your ergonomic considerations. You can adjust the height of the control panel to scan from a standing or sitting position, utilize the integrated footrest, and adjust the viewing angle of the monitor. For direct access to your frequently used functions, customize the touch screen buttons. Use the cable hangers to provide support for the transducer cables.

Available adjustments are dependent on the ultrasound system.

See also: For information about adjustments available for your ultrasound system, refer to Appendix A in this manual.

See also: Transducer Cable Management, page 3-17

System Startup


⚠ WARNING: Operating the ultrasound system in close proximity to other equipment can cause reciprocal interference. You should observe and ensure normal operation of the ultrasound system and other equipment.

The first step to operating the ultrasound system is to connect the system to a power supply outlet (mains outlet).

Plugging in the Ultrasound System

⚠ WARNING: Before connecting the ultrasound system to a power supply outlet (mains outlet), you must read and understand the Electrical Safety section of Chapter 2 in this manual.

The ultrasound system has a non-detachable power cord.

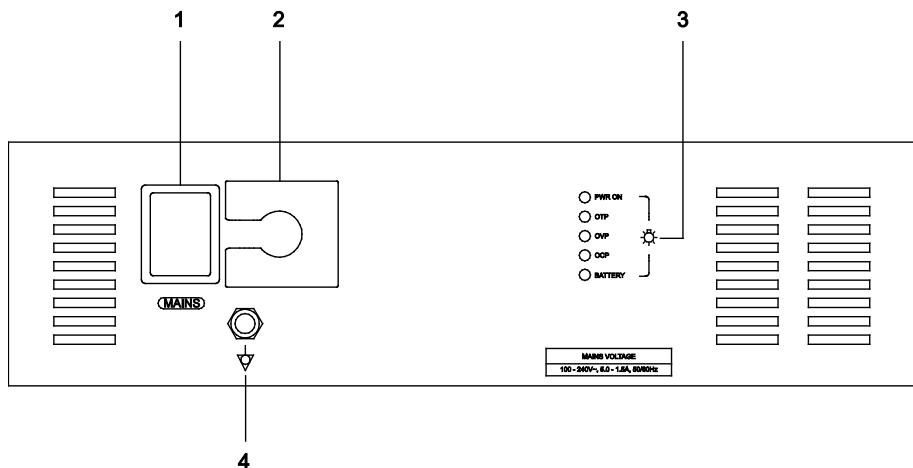

Example of power supply unit on the ultrasound system.

- 1 Mains switch
- 2 Connector cover for power cord
- 3 Status indicator lights for service diagnostics
- 4 Equipotential terminal

To plug in the ultrasound system:

- 200V~ to 240V~ systems: plug the power cord into a standard power supply outlet (mains outlet). For example, a "Schuko" outlet (CEE 7/7 standard).
- 100V~ to 120V~ systems used in the U.S.A.: plug the power cord into a hospital-grade power supply outlet (mains outlet).
- 100V~ to 120V~ systems outside the U.S.A.: plug the power cord into a standard power supply outlet (mains outlet).

Supplying Power to the Ultrasound System

The ultrasound system is powered on and off using the partial power on/off control located on the control panel.

Illumination of the Control	System Status
Blue	Powered on
Amber	Standby
Off (not lit)	Powered off

For optimal system performance, power off the ultrasound system every 24 hours.

To power on and off the ultrasound system:

 **Caution:** Wait approximately 20 seconds between powering the system off and then on again. This allows the system to complete its shutdown sequence.

Note: You cannot shut down and restart the ultrasound system when the battery support feature is activated.

1. To power on the system, press the partial power on/off control.
2. Visually check the on-screen displays and lighting as indicated in the daily checklist before using the ultrasound system.
See also: Refer to the daily checklist in Chapter 2 in this manual.
3. To power off the system, press the partial power on/off control.
 - If you configured the system to display a confirmation message, click **Shutdown**.
 - To cancel the shutdown, click **Close**.

To completely disconnect the ultrasound system from the power supply (supply mains):

1. Power off and unplug the ultrasound system from the power supply.
2. Ensure the circuit breaker on the ultrasound system (mains switch) is in the off position (O = OFF).

Activating System Standby

System standby decreases the time to power off and on the ultrasound system.

The ultrasound system preserves patient data when in system standby. To avoid the possibility of data loss, close the current patient study, activate 2D-mode, and wait until the system finishes writing data to external storage media before activating system standby.

Use the configuration settings to close the current study when system standby is activated.

System standby is not available during video playback or use of an application, for example, Stress Echo and Panoramic imaging. You cannot activate system standby with the battery support feature when the remaining level of battery charge is lower than 65%.

See also: For information about activating the battery support feature, refer to page 3-8.

To activate system standby:

1. Briefly press the partial power on/off control.
2. Click **Standby** and wait approximately 15 seconds for the standby feature to activate.
3. To power on the ultrasound system, press the partial power on/off control.

The system is ready for use in less than 15 seconds.

Activating QuikStart (Battery-Powered System Standby)

Activate system standby and power on the ultrasound system using the battery-powered QuikStart feature during a portable exam. When the system is disconnected from the power supply outlet, the available standby time is dependent on the installed battery pack.

- Battery packs only for QuikStart: approximately 50 minutes
- Battery packs for both QuikStart and battery support: approximately 180 minutes

The status of the battery charge is indicated by the lighting of the power on/off control.

Illumination of the Control	Battery Status
Steady amber	More than 60% charge remaining
Flashing amber (slow)	More than 50% charge remaining
Flashing amber (fast)	Battery power is low
Flashing amber (fast) with audible beep	Battery power is very low. If the ultrasound system is not plugged into the power supply outlet, the system shuts down within 5 minutes.

To activate QuikStart:

1. Briefly press the partial power on/off control.
2. Click **Standby**.
3. Wait approximately 15 seconds and then unplug the ultrasound system from the power supply outlet (mains outlet).
The ultrasound system activates the standby feature.
4. Plug the power cord into a power supply outlet (mains outlet).
5. To power on the ultrasound system, press the partial power on/off control.

Activating the Battery Support Feature

Use the battery support feature to scan a patient when the ultrasound system is unplugged from the power supply outlet and powered by the battery pack. The fully charged battery pack provides a maximum of 60 minutes of operation during 2D imaging without use of any documentation devices or accessories installed on the ultrasound system.

The battery support feature is not recommended for use during biopsy, stress echo, contrast, and 3D/4D volume imaging, or during intraoperative, transesophageal, or endocavity scanning.

Note: Perform the following tasks only when the ultrasound system is connected to the power supply outlet: installing system software or software updates, backing up or restoring configuration settings, and backing up or restoring patient data.

The ultrasound system displays icons and on-screen messages to indicate the level of battery pack charge, the recharging status of the battery, and any errors about the battery.

The number of bars on the on-screen icon indicates the level of battery pack charge. The color of the icon changes to yellow or red and the ultrasound system emits an audible beep when the battery reaches a low level of charge or is fully discharged. A power cord on the icon is visible when the ultrasound system is plugged into a power supply outlet and recharging the battery pack.



Example of the icon depicting a partially charged battery pack.



Example of the icon when the system is using power from the battery pack.

The capacity of a rechargeable battery degrades with use and over time. Expect an approximate 20% reduction in the original capacity after 200 full cycles of charging and discharging the battery. For more information on replacing the battery, contact your customer service representative.

To activate and use the battery support feature:

Note: You cannot shut down and restart the ultrasound system when the battery support feature is activated.

1. Plug the ultrasound system's power cord into a power supply outlet (mains outlet).
2. Power on the ultrasound system and wait approximately for 210 minutes to fully charge the battery pack.

Note: A battery pack can require approximately 210 minutes to fully charge at 25°C; however, the time required to fully charge a battery pack is dependent upon the existing battery pack charge and environmental factors, such as temperature and humidity. Expect an increase in the time to recharge a battery pack if the ultrasound system is not used daily or improperly shuts down. The capacity of a battery pack degrades over time and with use.

The following on-screen icon indicates a fully charged battery pack. The ultrasound system is plugged into a power supply outlet.



3. Unplug the ultrasound system from the power supply outlet (mains outlet) to provide power to the ultrasound system from the battery pack.

The following on-screen icon indicates a fully charged battery pack. The ultrasound system is unplugged from the power supply outlet and is using battery power.



4. Observe the status of the battery charge indicated by the battery icon on the image screen during system use.

The following on-screen icon indicates a critically low level of battery charge.



Troubleshooting On-Screen Messages

Follow the instructions included in the on-screen messages if the battery reaches a low battery level or exceeds the acceptable thermal limit.

Message	Displays when	Recommended Action
The system is out of battery. Connect the AC power.	The remaining level of battery charge is approximately 30%.	1. Observe the level of battery charge indicated by the on-screen icon during use of the ultrasound system.  2. Plug the power cord of the ultrasound system into the power supply outlet.
The system is out of battery. Connect the AC power. Or, the system will shut down in a few minutes.	The remaining level of battery charge is approximately 15%. The ultrasound system emits an audible beep.	1. Observe the level of battery charge indicated by the on-screen icon during use of the ultrasound system. A red icon blinks indicating the remaining battery charge is approaching a critically low level.  2. Immediately plug the power cord of the ultrasound system into the power supply outlet.
The system will shut down in a few seconds.	The remaining level of battery charge is approximately 0%. The ultrasound system emits an audible beep and can shut down within seconds.	1. Wait until the ultrasound system completely shuts down.  2. Plug the power cord of the ultrasound system into the power supply outlet to recharge the battery pack.
The battery is overheated. Wait until the battery temperature reaches a normal operation level with the AC power connected.	The battery temperature exceeds the normal operation temperature limit.	Note: Do not unplug the power cord of the ultrasound system from the power supply outlet. 1. Stop using the ultrasound system when the battery pack is overheated and the ultrasound system displays the following on-screen icon.  2. Wait until the ultrasound system removes the exclamation mark (!) from the on-screen icon and then continue to use the ultrasound system.

Recalibrating the Battery Pack

Recalibrate the battery pack every 180 days to maintain optimal battery performance. The ultrasound system is not available for use during the recalibration process. Battery recalibration requires several hours. Prior to recalibrating the battery pack, consult with your customer service representative.

To recalibrate the battery pack:

1. Plug the ultrasound system's power cord into a power supply outlet (mains outlet).
2. Power on the ultrasound system and wait approximately for 210 minutes to fully charge the battery pack.

Note: A battery pack can require approximately 210 minutes to fully charge at 25°C; however, the time required to fully charge a battery pack is dependent upon the existing battery pack charge and environmental factors, such as temperature and humidity. Expect an increase in the time to recharge a battery pack if the ultrasound system is not used daily or improperly shuts down. The capacity of a battery pack degrades over time and with use.

The following on-screen icon indicates a fully charged battery pack.



3. Do not use the ultrasound system for an additional two hours and then unplug the ultrasound system from the power supply outlet.
4. Wait until the battery pack is completely discharged and the ultrasound system shuts down.

The following on-screen icon indicates a completely discharged battery pack.




5. Do not use the ultrasound system for an additional two hours.
6. Repeat steps 1 and 2 to fully recharge the depleted battery pack.

Connecting and Disconnecting Transducers

⚠ Caution: You must freeze the system before connecting or disconnecting a transducer. If a transducer is disconnected before the image is frozen, the system will display an error message, and it will be necessary to reset the system before continued use.


⚠ Caution: To avoid damaging the transducer cable, do not pull on the cable to disconnect a transducer.

ACUSON Juniper System

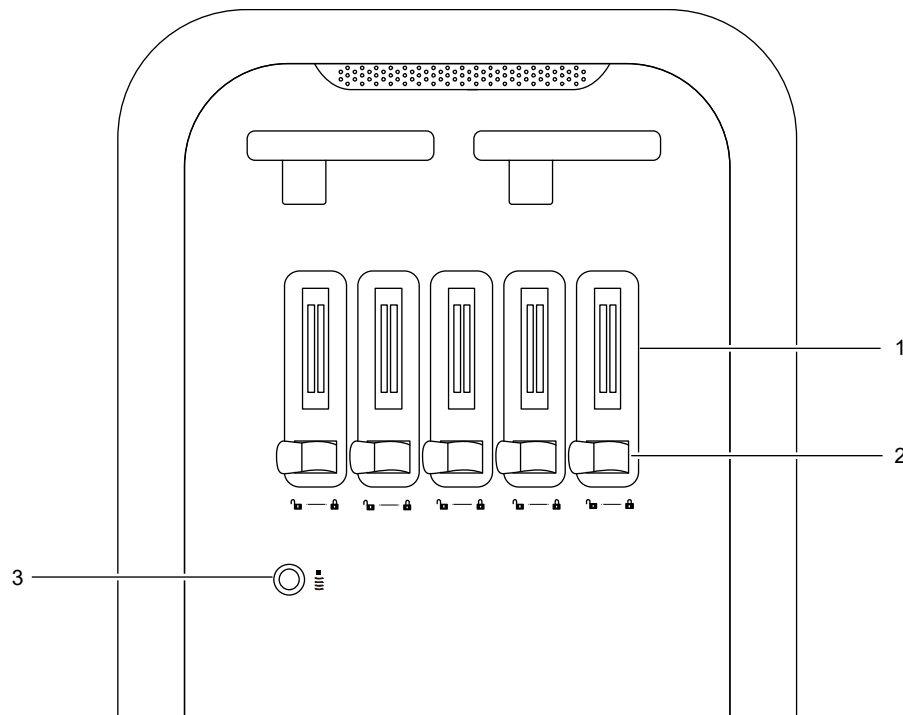

Example of transducer ports.

- 1 Transducer ports for array transducers
- 2 Transducer port locking lever
- 3 Transducer port for continuous wave (pencil) transducers

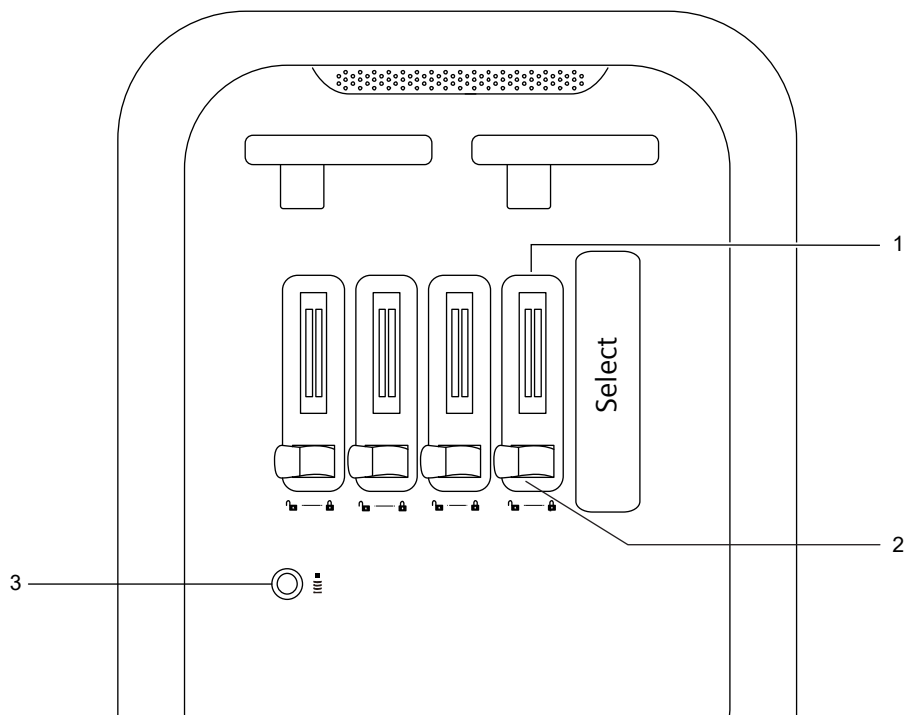ACUSON Juniper Select System

Example of transducer ports.

- 1 Transducer ports for array transducers
- 2 Transducer port locking lever
- 3 Transducer port for continuous wave (pencil) transducers

Array Transducers

When array transducer connectors are connected or disconnected from the ultrasound system, there may be some resistance due to the shielding material inside the connectors. This is normal for these transducers.

To connect an array transducer:

1. Make sure the locking lever below the transducer port is in the unlocked position.



2. Insert the transducer connector into the port, taking care to align the connector for a smooth fit.
3. Push the locking lever to the locked position.



The system detects the transducer and displays the transducer name on the image screen and touch screen.

To disconnect an array transducer:

1. Push the locking lever to the unlocked position.



2. Firmly grasp the transducer connector and carefully pull it straight out from the transducer port.

See also: For information about cleaning, disinfecting, and storing transducers, refer to Chapter 6 in this manual.

Continuous Wave Transducers

Use the transducer port for continuous wave (pencil) transducers.


To connect a continuous wave transducer:

1. Align the connector key until the connector fits smoothly into the receptacle.
2. Insert the connector into the transducer port until the connector locks in position.

To disconnect a continuous wave transducer:

- Pull on the connector housing ring to disengage the locking mechanism. Continue pulling on the ring to remove the connector from the transducer port.

See also: For information about cleaning, disinfecting, and storing transducers, refer to Chapter 6 in this manual.

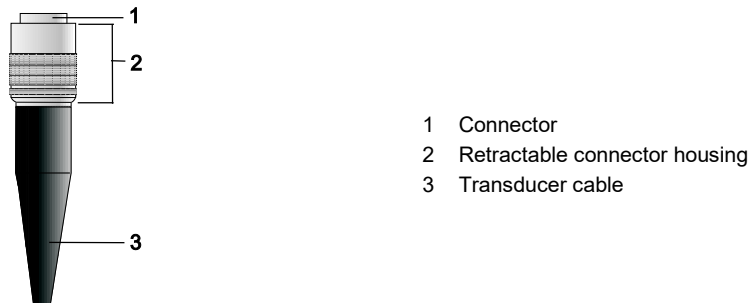


Example of the connector on a continuous wave transducer.

Protective Transducer Holders and Gel Holders

⚠ Caution: Transducer holders have variable sizes both in depth and diameter. To avoid transducer damage, you must use the holder and/or insert provided for specialty transducers with small or large diameter handles.

After connecting a transducer to the ultrasound system, place the transducer in the protective holder attached to the control panel. Holders are also available for standard-size and large-size gel bottles.

The transducers holders are interchangeable and replaceable.

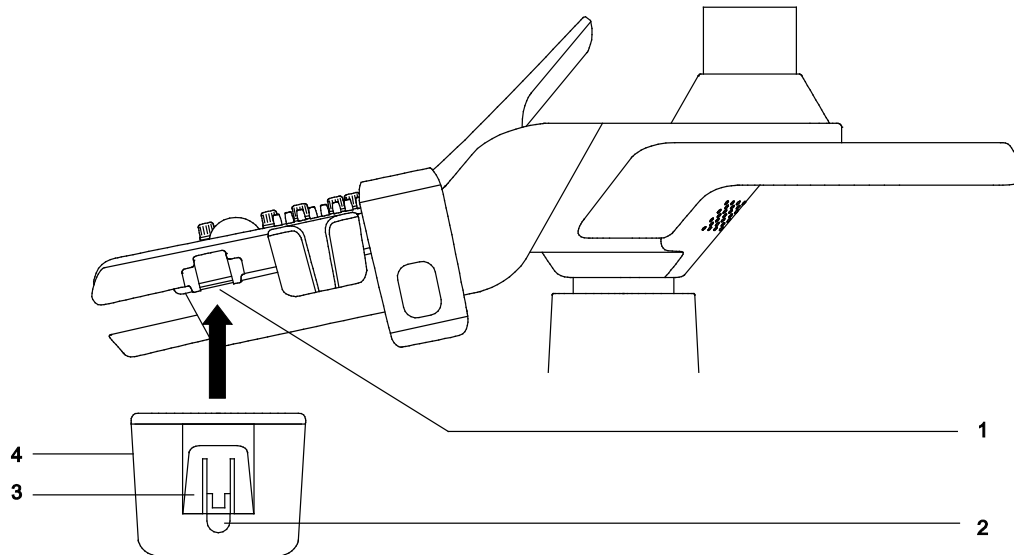

Example of transducer holders.

- 1 Point of attachment
- 2 Tab
- 3 Support
- 4 Holder

Note: The gel holder for the large bottle uses the same point of attachment as the transducer holders.

To install a holder:

1. Press the tab against the holder.
2. Align the support on the holder with the attachment slot on the system and slide the holder upward until it snaps into place.

To remove a holder:

1. Locate the tab underneath the holder.
2. Press the tab against the holder and pull the holder downward.

Transducer Cable Management

 **WARNING:** To avoid injury from tripping over a transducer cable, use the cable hooks on the ultrasound system to keep the transducer cable from contact with the floor.

After you connect and secure a transducer, drape the transducer cable through one of the cable hooks located on the ultrasound system. These hooks provide support for the transducer cables, keep cables off the floor, and help to prevent tangling of the cables when more than one transducer is connected to the system.


Gel Warmer

⚠ WARNING: To avoid electrical shock and damage to the control panel resulting from ingress of liquid, place the gel bottle and gel warmer on the side of the ultrasound system closest to the patient.

⚠ Caution: Do not attempt to repair or alter the gel warmer. Contact your local Siemens Healthineers service representative immediately if your gel warmer appears to be damaged or malfunctions in any way.


See also: For information about cleaning and disinfecting the gel warmer, refer to Chapter 2 in this manual.

When the ultrasound system and the gel warmer are powered on, the gel warmer continuously heats the bottle of gel inserted into the gel warmer.

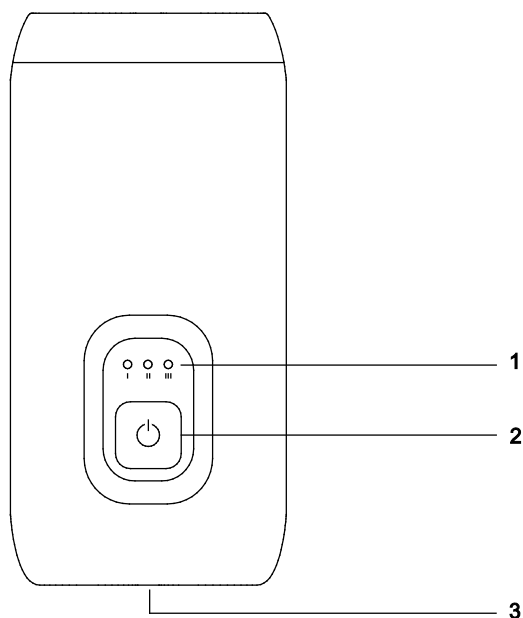

Example of the gel warmer, front view.

- 1 Temperature indicator
- 2 Control for power and temperature
- 3 Cap for the inverted gel bottle

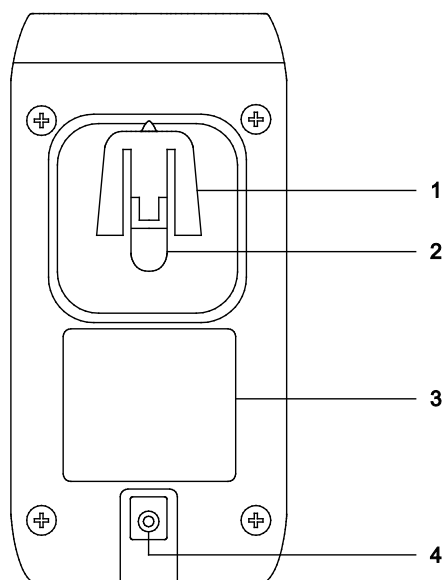

Example of the gel warmer, back view.

- 1 Support
- 2 Tab
- 3 Label
- 4 Power cable connector

To attach the gel warmer to the ultrasound system:

1. Remove the transducer holder on the side of the ultrasound system closest to the patient and farthest from the front of the system.
 - a. Reach under the holder to locate the tab on the holder. The tab extends below the point of attachment to the ultrasound system.
 - b. Squeeze the tab towards the holder and pull the holder downward.
2. Attach the gel warmer to the ultrasound system.
 - a. Attach one end of the power cable accompanying the gel warmer to the connector on the back of the gel warmer and then attach the other end of the power cable to the connector on the ultrasound system.
 - b. Align the support on the gel warmer directly below the point of attachment on the ultrasound system and firmly push upwards until the gel warmer snaps into place.

To use the gel warmer:

1. Press the control on the front of the gel warmer to power on the gel warmer, if necessary. The gel warmer power control illuminates when the gel warmer is powered on.

Note: A flashing light indicates there is a problem with the device. Check the cable connector on the gel warmer and on the ultrasound system. If you continue to experience a problem with the device, contact your Siemens Healthineers service representative.

2. Press the control on the front of the gel warmer to adjust the temperature of the gel warmer.

Selection	Setting	Temperature
I	Lowest	31°C
II	Medium	34°C
III	Highest	37°C


3. Invert the bottle of gel, and insert the bottle into the gel warmer.

Installing the Storage Container

You can attach or remove the basket or storage bin on the system.

To install the basket on the ultrasound system:

1. Attach the basket on the rear handle of the system.

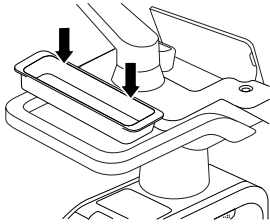


Example of attaching the basket on the rear handle.

2. To remove the basket from the system, gently lift the basket vertically from the rear handle.

To install the storage bin on the ultrasound system:

1. Attach the storage bin on either the right side or left side of the system.
An audible click indicates the storage bin is properly attached to the system.

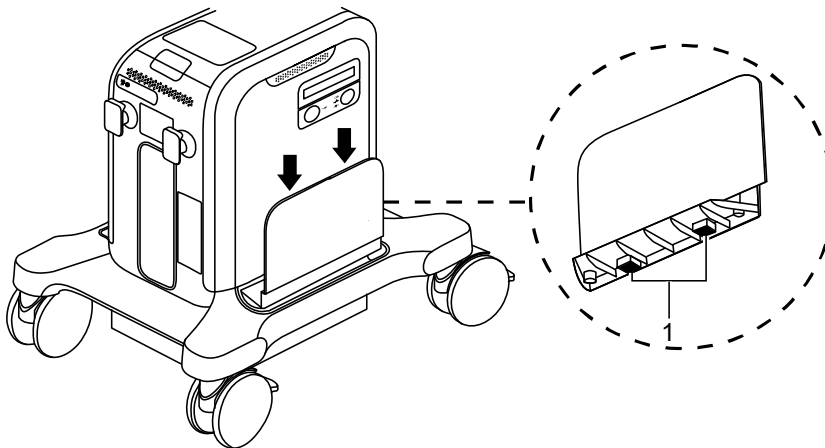

Example of attaching the storage bin on the side of the system.

- 1 Magnets at the bottom of the storage bin for attachment to the system
2. To remove the storage bin from the system, gently lift the storage bin vertically from the system.

Adjusting the Audio Volume

You can adjust the Doppler volume and adjust the volume of audible tones (beeps) that confirm key presses and store functions.

Use the configuration settings to adjust the audio volume and enable or disable confirmation beeps.

To adjust the volume settings:

Prerequisite: Assign adjustment of the volume to the touch screen controls in the configuration settings.


1. Rotate **[Volume]**.
2. To mute the volume, tap **Mute** or rotate **[Volume]** to turn off the volume.

Adjusting the Screen Brightness

You can increase or decrease the brightness of the touch screen and adjust settings on the monitor for the image screen. Make adjustments to the brightness and contrast of the monitor before adjusting the print quality of installed documentation devices.

Note: Factory-defined imaging presets were created using default settings for the brightness of the monitor. Adjusting the brightness on the monitor may affect the image optimization intended by the factory-defined image presets.

Note: The language setting on the monitor menu is separate from the language setting of the ultrasound system.

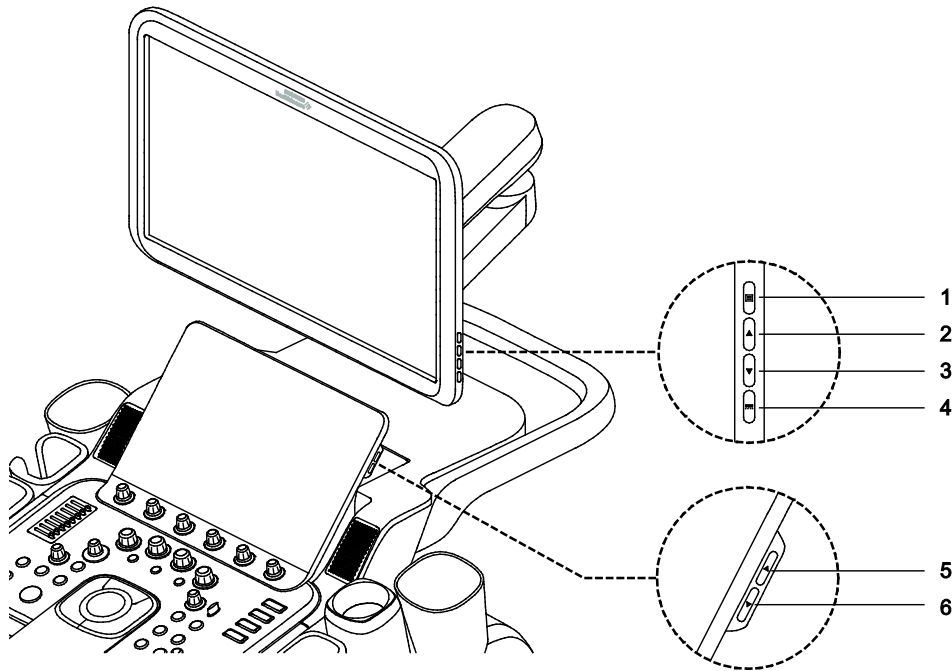

Example of monitor and touch screen controls.

- 1  accesses the menu for the monitor
- 2  adjusts settings on the monitor
- 3  adjusts settings on the monitor
- 4  accesses brightness settings for the lamp brightness
- 5  increases brightness of the touch screen
- 6  decreases brightness of the touch screen

To adjust the brightness of the monitor:

1. Press the monitor menu control.



2. Press the up or down controls to adjust the setting.



Note: To restore the default contrast setting, simultaneously press the up and down controls.

3. To exit the menu, press the monitor menu control.

To adjust the brightness of the touch screen:

- Press the up or down controls to adjust the setting.



Input/Output Connections


⚠ WARNING: Accessory equipment connected to the analog and digital interfaces must be certified according to the respective EN or IEC standards (for example, EN 60950 or IEC 60950 for data processing equipment and EN 60601-1 or IEC 60601-1 for medical equipment). Anyone who connects additional equipment to any of the signal input or signal output ports configures a medical system and is therefore responsible that the system complies with the requirements of the system standards EN 60601-1 or IEC 60601-1. Siemens Healthineers can only guarantee the performance and safety of the devices listed in the *Instructions for Use*. If in doubt, consult the Siemens Healthineers service department or your local Siemens Healthineers representative.

See also: For additional information about accessories and options, refer to Appendix A in this manual.

⚠ WARNING: To ensure proper grounding and leakage current levels, it is the policy of Siemens Healthineers to have an authorized Siemens Healthineers representative or approved third party perform all on-board connections of documentation and storage devices to the ultrasound system.

⚠ Caution: To prevent damage to the ultrasound system and peripheral devices, always power off the system before connecting or disconnecting peripheral devices.

Network, audio, and video connections are located on the input/output panel.

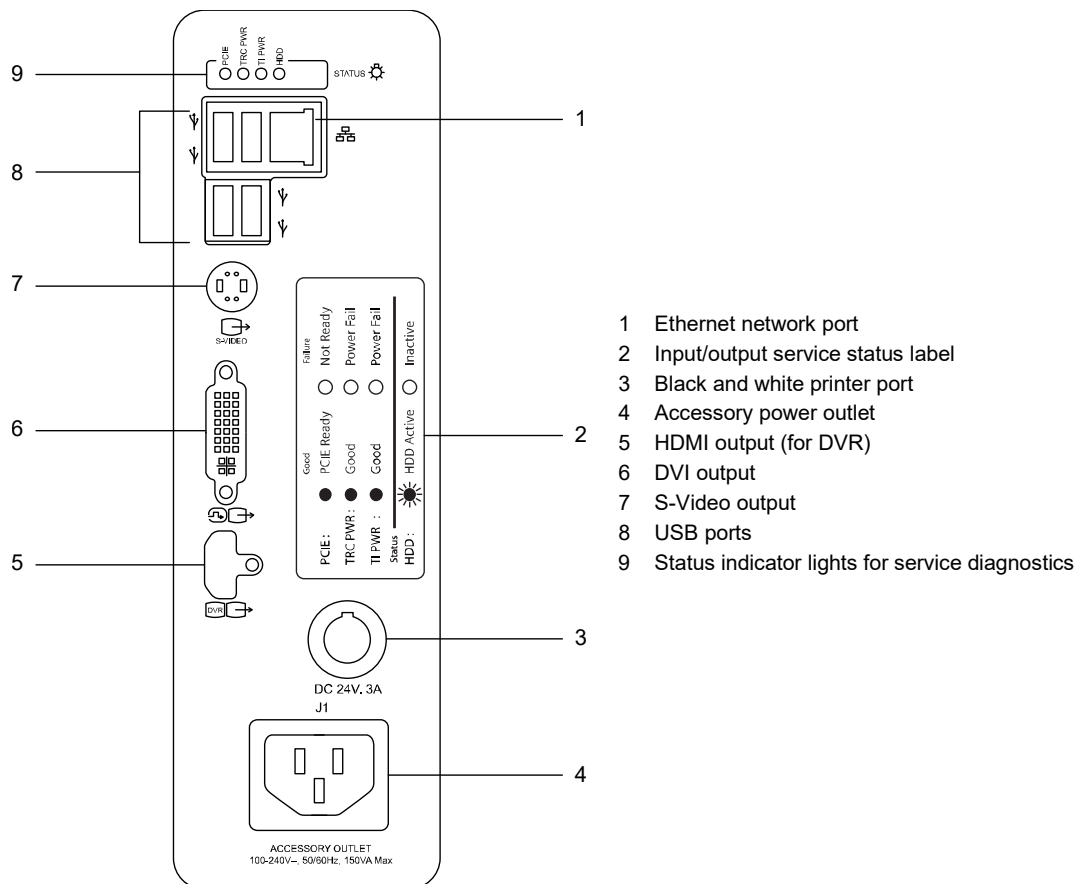

Example of input/output panel.

Connecting Peripheral Equipment and Accessories

 **WARNING:** Equipment connected to the ultrasound system and in the patient environment must be powered from a medically-isolated power supply outlet (mains outlet) or must be a medically-isolated device. Equipment powered from a non-isolated source can result in chassis leakage currents exceeding safe levels. Chassis leakage current created by an accessory or device connected to a non-isolated outlet may add to the chassis leakage current of the ultrasound system.

 **WARNING:** Non-medical grade report printers cannot be used within a patient environment.

 **WARNING:** During use of a non-medical grade report printer or when a non-medical grade report printer is connected to the ultrasound system, the ultrasound system cannot in any way be in contact with a patient.

Peripheral devices must be installed on the ultrasound system by an authorized Siemens Healthineers representative or by a Siemens Healthineers approved third party. Any use of other devices with the system will be at the user's risk and may void the system warranty.

Connect peripheral devices and accessories to the USB ports located on the ultrasound system. For example, you can connect a medical grade printer, a medical grade DVR, or a USB storage device to the ultrasound system.

Use the configuration settings for the ultrasound system to add a printer to the list of devices recognized by the ultrasound system.

To fulfill EN 60601-1 or IEC 60601-1 (Medical Electrical Equipment, Part 1: General Requirements for Basic Safety and Essential Performance) requirements, connection of peripheral equipment to your ultrasound system must adhere to one of the following conditions:

- The peripheral equipment itself is a medical device approved according to EN 60601-1 or IEC 60601-1; or
- Non-medical peripheral equipment approved according to any other EN or IEC standard (EN XXXXX or IEC XXXXX, for example, equipment complying with EN 60950 or IEC 60950) must use the following setup for connection:



Equipotential terminal located on the power supply unit (PSU) of the ultrasound system.

- Connect the ultrasound system to an independent protective ground terminal, with a ground wire connection to the ultrasound system's equipotential connector. Ensure that the protective ground wire is connected to a qualified protective ground connection independent of the existing system's ground connection (via the power cord).
- The peripheral equipment is located at least 1.5 meters (1.8 meters [6 feet] in Canada and the U.S.A.) outside of the patient environment. A patient environment is defined as the area in which medical examination, monitoring, or treatment of the patient takes place.
- The peripheral equipment is connected to a mains outlet outside the patient environment, but still within the same room as the ultrasound system.

For additional information and other possible combinations, refer to the Medical Electrical Equipment Standard EN 60601-1 or IEC 60601-1, Annex I.

Note: The above information is based on current EN 60601-1 or IEC 60601-1 standards. If your country's regulatory standards for medical equipment do not correspond to EN 60601-1 or IEC 60601-1, your local requirements may be different.


To connect a peripheral device to the ultrasound system:

1. Insert the connector from the peripheral device into the appropriate port on the input/output panel.


Note: If a printer is not recognized by your ultrasound system, use the configuration settings for the ultrasound system to add a printer to the list of recognized devices.

2. Observe the peripheral equipment connections and patient environment shown in the following diagrams to support the normal operation of the ultrasound system.

Note: During normal operation, the ultrasound system is designed to display ultrasound images without noise, artifacts, or distortion that cannot be attributed to physiological effects.

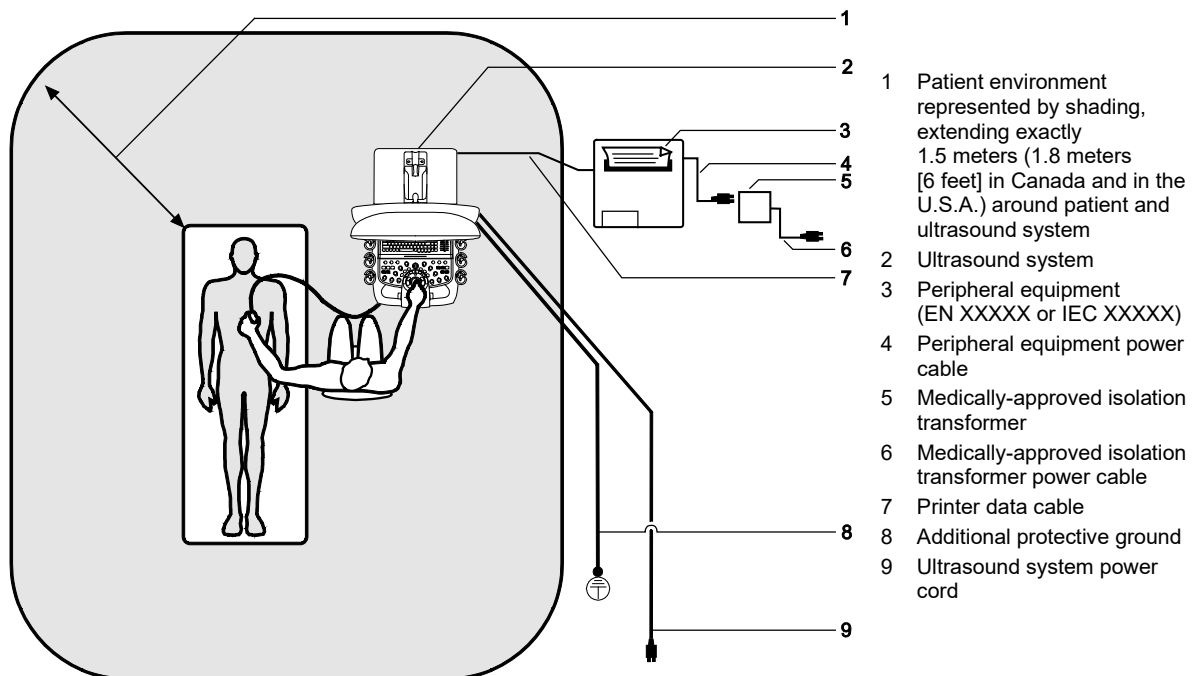

Example of a peripheral equipment connection and patient environment.

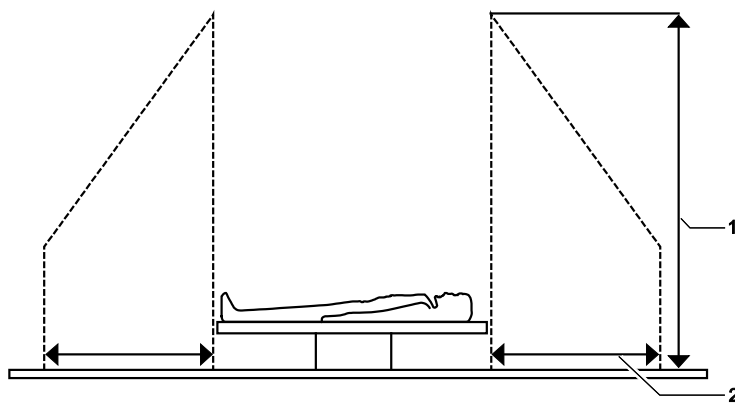

Example of a peripheral equipment connection and patient environment.


- 1 2.5 meters (typical)
- 2 1.5 meters (1.8 meters [6 feet] in Canada and the U.S.A.)

Installing the Protective Overlay

The protective overlay is an accessory to avoid possible contamination to the control panel, rotary controls, and touch screen.


To install the protective overlay on the control panel:

1. Cover the control panel with the protective overlay.

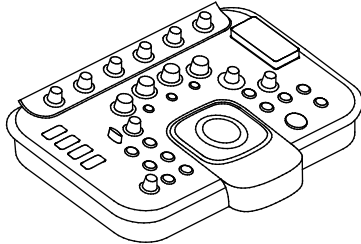


Example of overlay for the control panel.

2. Place the following protective overlays over the controls.

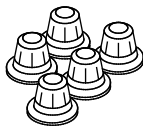

Example of overlays for rotary controls.

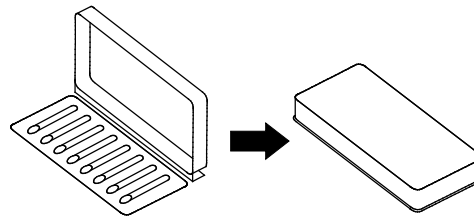

Example of overlay for the depth gain compensation controls.

3. Cover the touch screen with the protective overlay.

Connecting the Ultrasound System to a Network

You can connect the ultrasound system to a local area network (LAN) through a wired or a wireless connection. When both wired and wireless connections are configured and available, the ultrasound system assigns priority to the wired connection.

If you are in the process of saving a study to a network location when you disconnect the system from a LAN or change the connection to a LAN, the ultrasound system pauses the process of saving the study and resumes saving the study when you reconnect the system to the LAN.

Wired Connection

To connect the ultrasound system to a network through a wired connection:

Prerequisite: The network administrator has configured the required wired connection on the ultrasound system.

1. Connect the network cable to the Ethernet port located on the input/output panel.
2. Connect the network cable to the hospital or clinic network port.

If the ultrasound system is connected to a wired network, the ultrasound system disconnects from the wireless network. The system then connects to the LAN through the wired connection and displays an icon in the system status bar indicating a wired connection.



Example of the icon indicating a wired connection.

Wireless Network Connections

The ultrasound system can send data, such as studies, images, and clips, to a network location over a wireless network.

The ultrasound system displays an icon indicating the status and signal strength of the wireless network in the system status bar.

Icon	Description
	Indicates the signal strength of the wireless network.
	Indicates the wireless network is not transmitting a signal. Use the configuration settings to enable connections to a wireless network.

Use of Wireless Adapters

 **WARNING:** During use of a non-medical grade wireless adapter or when a non-medical grade wireless adapter is connected to the ultrasound system, the ultrasound system cannot in any way be in contact with a patient.

 **WARNING:** To prevent the loss of data over a wireless network and to avoid damaging the ultrasound system, never use a wireless adapter not recommended by Siemens Healthineers.

For information on the transmission and reception frequencies used by the wireless adapter, refer to the manufacturer's operating instructions included with the wireless adapter.

You must abide by the regulatory requirements of your country to use a wireless adapter. Refer to the manufacturer's operating instructions for information regarding certified use of the wireless adapter in your country.

For information on the compatible wireless adapters, contact your customer service representative.

Connecting to Remote Services

(Available only when remote services is enabled on your ultrasound system)

You can connect to remote services through a wired or wireless network.

To begin the remote access session:

1. Use the phone number provided by your Siemens Healthineers representative.
2. Follow the instructions provided by your service representative, including how to disconnect from the remote service session.

Note: If you do not disconnect from the session, the session will end after 30 minutes of inactivity.

Moving the System

-  **WARNING:** Preparations before moving the ultrasound system are important to minimize potential damage to sensitive components and to avoid safety hazards. Review the moving instructions before moving the system.
-  **WARNING:** Use caution when moving or positioning the ultrasound system. Exerting excessive force can cause the system to lose balance, creating a risk for human injury and damage to the system.
-  **Caution:** Do not park the ultrasound system or leave the system unattended on a slope. Even when the rear wheels are locked, the system may slide down a ramp.
-  **Caution:** Do not lean on the monitor, touch screen, storage bin, or keyboard. Subjecting these components to heavy loads or excessive force can damage the ultrasound system.

The ultrasound system is a mobile unit. Before moving the system to another location, you must prepare for the move by powering off and securing the system.


Securing the Wheel Locks

The ultrasound system has locking wheels. Use your foot to lock or unlock each wheel.

Wheels with a directional swivel lock are not available on all ultrasound system models.

See also: For information about each system model, refer to Appendix A in this manual.

Wheels without Swivel Locks

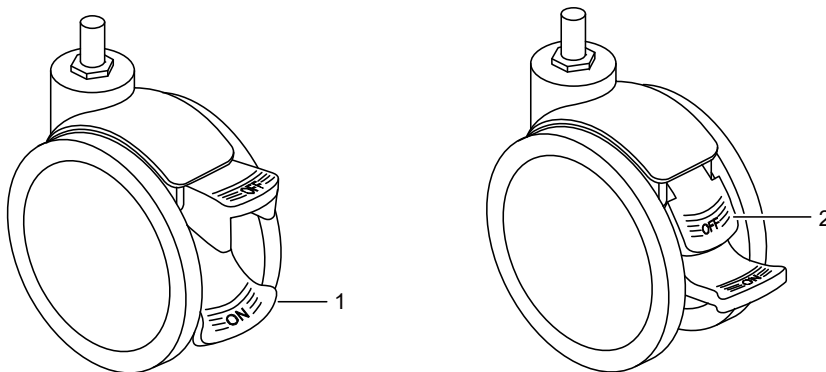


Example of a wheel with a locking lever and unlocking lever.

- 1 Locked wheel
The **ON** lever is in a downward position.
- 2 Unlocked wheel
The **OFF** lever is in a downward position.

To lock and unlock the wheels:

- To prevent the wheels from rolling, push downward on the **ON** levers.
- To unlock the wheels, push downward on the **OFF** levers.

Wheels with Swivel Locks

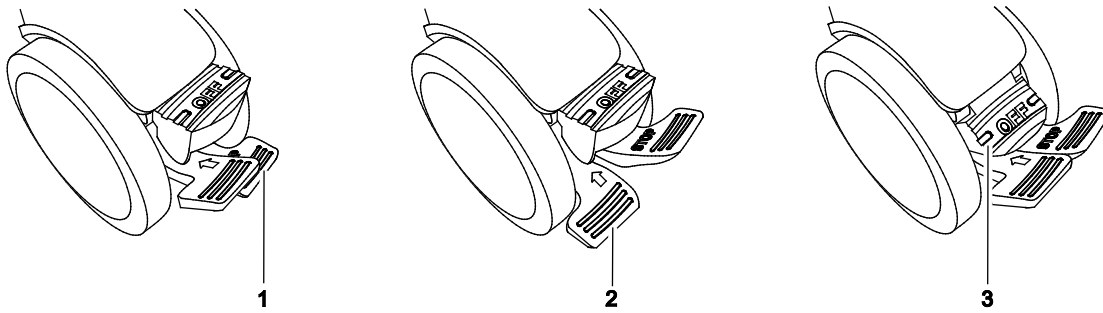

Example of a wheel with a locking lever, swivel lock lever, and unlocking lever.

- 1 Locked wheel
The **STOP** lever is in a downward position.
- 2 Swivel locked wheel
The arrow lever is in a downward position.
- 3 Unlocked wheel
The **OFF** lever is in a downward position.

To lock and unlock the wheels:

- To prevent the wheels from rolling and swiveling, push downward on the **STOP** levers.
- To prevent the wheels from swiveling, push downward on the arrow levers.
- To unlock the wheels, push downward on the **OFF** levers.

Before the Move


-  **WARNING:** To avoid damage to the monitor and the potential for personal injury to the user, ensure the monitor does not swivel during transport. Ensure the monitor is in the vertical, locked position prior to moving the ultrasound system according to the directions provided in this section.
-  **WARNING:** Carefully position the monitor to avoid risk for injury to the user. A potential pinch point between the adjustable arm and the touch screen exists when swiveling the adjustable arm to the 180° position. A potential pinch point also exists between the monitor and the adjustable arm when rotating the monitor.
-  **Caution:** To prevent damage to the system during a move, retract or close any doors and trays on documentation devices or components on the ultrasound system. Ensure components do not protrude from the system.

Always pull on the plug when disconnecting the ultrasound system from the power supply outlet (mains outlet) and disconnecting devices from the system. Never pull on cords or cables.

ACUSON Juniper System

To prepare the ultrasound system for the move:

1. Power off the ultrasound system.
2. Unplug the power cord from the power supply outlet (mains outlet).
3. Secure the power cord to avoid rolling the system wheels over the cord.
4. Disconnect all devices from the system, including the network cable, off-board peripheral devices, system accessories, and transducers.
5. Place each clean transducer in a protective carrying case for safe transport.
6. Close any doors and trays on documentation devices or components on the ultrasound system and retract the keyboard, if necessary.
7. Transport gel and disk media separately.
8. Position and lock the monitor for transport:
 - a. Align the monitor to the front, center of the ultrasound system.
 - b. Rotate the transport lock into the locked position. Ensure the locking pin engages with the hole on the adjustable arm.

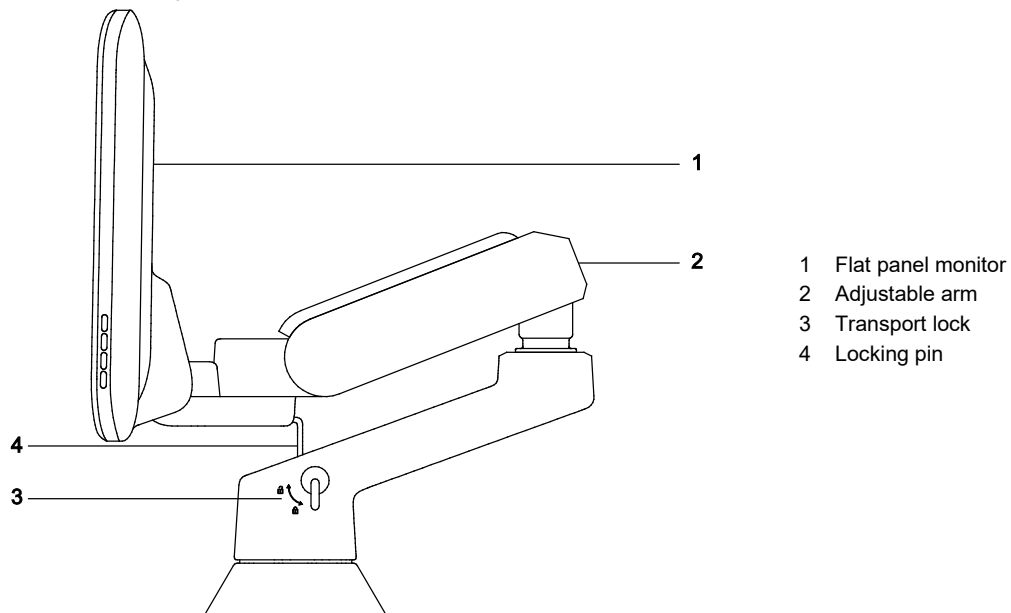


Example of monitor in vertical position for transport.

9. Unlock the front and rear wheels.

ACUSON Juniper Select System

To prepare the ultrasound system for the move:

1. Power off the ultrasound system.
2. Unplug the power cord from the power supply outlet (mains outlet).
3. Secure the power cord to avoid rolling the system wheels over the cord.
4. Disconnect all devices from the system, including the network cable, off-board peripheral devices, system accessories, and transducers.
5. Place each clean transducer in a protective carrying case for safe transport.
6. Close any doors and trays on documentation devices or components on the ultrasound system and retract the keyboard, if necessary.
7. Transport gel and disk media separately.
8. Align the monitor to the front, center of the ultrasound system.

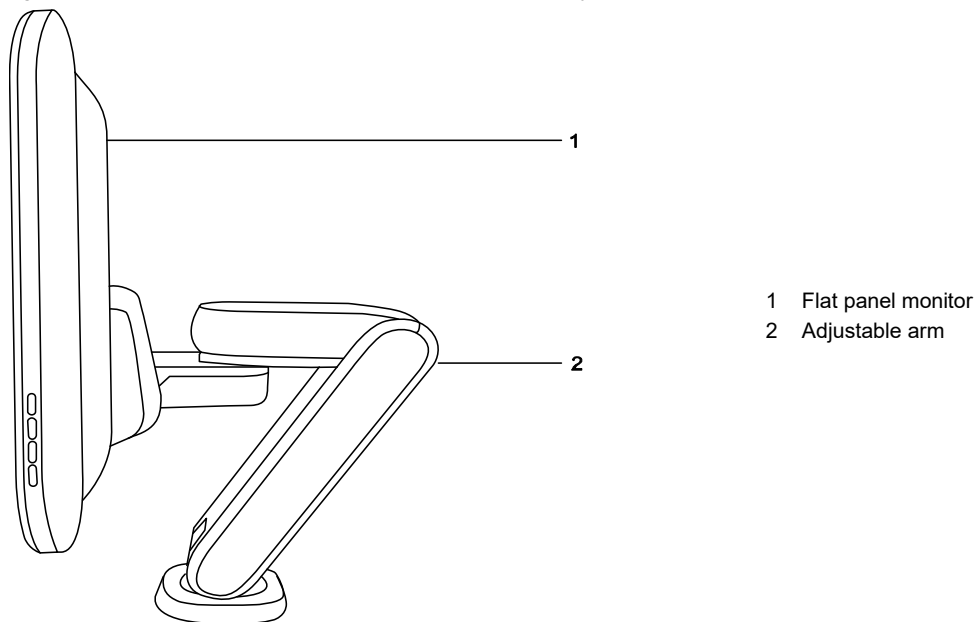

Example of monitor in vertical position for transport.

9. Unlock the front and rear wheels.

During the Move

-  **Caution:** When moving the ultrasound system, protect the system from environmental changes including moisture, winds, dirt and dust, and extreme heat or cold exposure.
-  **Caution:** Avoid moving the ultrasound system on outside surfaces with loose dirt, contaminants, or standing liquids.
-  **Caution:** Care should be taken to minimize shock and vibration of the ultrasound system. Avoid uneven surfaces that contain an abrupt height change or jarring surface irregularities.
-  **Caution:** For systems installed with an on-board documentation device located on the rear shelf of the ultrasound system, avoid damage to cables and connectors protruding from the documentation device, particularly when moving the system around corners or through doorways.
-  **Caution:** Do not move the ultrasound system by pushing on the monitor, touch screen, storage bin, or cable hooks. Pushing on these components can cause loss of control and damage the ultrasound system. Use the front handle to move or reposition the system.

You can move the ultrasound system from room to room within a facility and reposition the system during an examination. Be careful on inclines and uneven surfaces. The ultrasound system can be moved across pavement and other hardened parking lot surfaces.

Note: You must lock the wheels of the ultrasound system when transporting the system by vehicle. The ultrasound system must be sufficiently anchored to the vehicle floor or walls so the system to avoid shifting or movement during transport.

After the Move

-  **Caution:** Obstructed fans can cause potential ultrasound system overheating, system performance degradation, or failure. Do not allow fabric, such as linens or bedding, or surfaces to impede ventilation around the system. Ensure the system has proper ventilation during operation.
-  **Caution:** Wheel locks are most effective on a level surface. Never park the ultrasound system on an incline greater than five degrees.

To set up the ultrasound system after a move:

1. Position the ultrasound system: Make sure the system is not placed against walls or fabrics that obstruct perimeter air flow to the system cooling fans.
2. Lock the front and rear wheels.
3. Plug the power cord into a hospital-grade or local equivalent power supply outlet (mains outlet).
4. Reconnect devices to the system, including the network cable, off-board peripheral devices, system accessories, and transducers.
5. Power on the ultrasound system.

Shipping the Ultrasound System

You must prepare the system for shipment over long distances or rough terrain.

To prepare the ultrasound system for shipment:

1. Repack the ultrasound system in the factory packaging and crate.
2. Load the system into a vehicle using a lift gate.
3. Secure the system with cargo straps to prevent lateral movement of the system during transport.
4. Provide shock cushions beneath the system to prevent sudden jarring of the system during transport.

Installing Updates

(Available only when remote services is enabled on your ultrasound system)

The remote update handling feature sends updates to your ultrasound system over a secure network connection. You should install an update within one week after the update is downloaded. However, you cannot postpone a mandatory security update.

Note: Installing a software update requires a permission level of administrator.

To install or defer updates:

Prerequisite: Ensure the system is connected to a wired or wireless network.

1. Power on the ultrasound system.
The system displays the update installation window.
2. Click **Description** to view a description of the update.
3. Select the check box to confirm you have read the description.
4. Click **Install**.

The update description provides the required installation time.

If multiple updates are ready for installation, the system installs the first update and then prompts you to continue installation of the next update.

5. Click **OK**.

Note: Do not power off the ultrasound system while installation is in progress.

The system may restart multiple times during installation. The system displays a message when installation is complete.

6. To postpone installation of the update, click **Install Later**.
7. To postpone installation of all updates, click **Install All Later**.

Installing System Software

You can install new system software or reinstall existing system software with or without keeping your user-defined settings and patient data.

Note: Siemens Healthineers recommends that you retain user-defined settings and patient data. Whether you retain or remove user-defined settings and patient data, you should make a backup copy of system configuration settings and user-defined exams before installing software.

See also: For information about storing and retrieving system configurations and user-defined exams, refer to Chapter 1 in the System Reference.

Use the configuration settings to verify the system software version.

To install the system software:

Prerequisite: If your system is not equipped with the retractable keyboard, connect a keyboard to the USB port. The touch screen keyboard is not available during software installation.

1. Connect the installation media to the USB port.
2. Briefly press the partial power on/off control.
3. Keep pressing and releasing **F11** on the keyboard until the system displays the window for BIOS password input.
4. Enter the BIOS password provided by your Siemens Healthineers representative and then press enter on the keyboard.

The system displays the boot menu to select the installation media.

5. Select the installation media inserted and then press enter on the keyboard.
The system displays the software installation window.
6. Specify the required settings for the software installation:
 - **Language** specifies the language for the operating system.
 - **Retain patient data** (available only for reinstallation) saves and restores the patient data even after installation.
 - **Retain user data** (available only for reinstallation) saves and restores the user-defined configuration settings even after installation.

7. Click **Apply** and then click **OK** when the system displays a confirmation message.
The system restarts when installation is complete.

Note: Do not power off the ultrasound system while the installation is in progress.

8. Remove the installation media from the USB port.

Installing a Software-Based Option

To install a software-based option:

1. Tap the configuration icon.



2. Tap **Service**.
3. Click **Service** to access the settings in the service software and then click **OK** to display the user-configurable settings.
4. Click **Configuration**.
5. Click **Licensing**.
6. Connect the installation media to the USB port.
7. Click **Select New License File**.
8. Click the license file in **Browse** or **History**.
9. Click **Save**.

The software-based options are displayed in the **Available Features** section.

10. Click **Home**.
The system displays a restart dialog box.
11. Click **OK**.
The system restarts and confirms the installation of the options.

Accessing Help

The User and Reference Manuals for your ultrasound system are available from the Siemens Healthineers Document Library at doclib.siemens-healthineers.com.

To locate your user manual in the online document library:

1. Tap the help icon.



The system displays the single-page instruction with a QR code.

2. Scan the QR code on the ultrasound system with your mobile device.
3. Open the Siemens Healthineers Document Library on your mobile device.
Note: Register for an online account, if necessary.
4. Filter or browse the library to locate the user and reference manual for your ultrasound system.
5. Open the user manual and search for information.

4 Examination Fundamentals

Preventing Unauthorized Access.....	3
Logging In to the Ultrasound System.....	3
Registering a Patient and Beginning a Study	5
Preregistering a Patient	5
Registering a New Patient and Beginning a Study	5
Selecting an Existing Patient and Beginning a Study	6
Registering a Patient from the Worklist Server.....	7
Customizing the Patient Registration Form	8
Editing Patient Registration Information	8
Beginning a Study without Registering the Patient.....	8
Restarting an Existing Study.....	9
Selecting or Customizing the Exam Type	10
Customizing Exams	11
Activating and Deactivating Operating Modes.....	12
Workflow Protocol	13
Optimizing Images.....	14
Adjusting the Settings for General Imaging Parameters.....	14
2D Imaging	15
Optimizing 2D Images	15
Color Imaging	18
Acoustic Streaming with Cysts	19
Optimizing Color Images	19
Pulsed Wave Doppler.....	22
Optimizing the Doppler Spectrum.....	23
M-mode	25
Optimizing the M-mode Sweep	25
Examples of Imaging Settings.....	26
Settings for All Modes.....	26
2D Settings	26
Color Settings	26
Pulsed Wave Doppler Settings	27
M-mode Settings	27
General Contrast Agent Imaging Settings	27

CINE 28

Capturing Images and Clips 29

Reviewing Images and Clips 29

Printing Images..... 29

Video Recording a Study 30

Completing an Examination 30

Preventing Unauthorized Access

The security package on the ultrasound system protects patient information from unauthorized access by requiring a user account and password to access the ultrasound system. It is the responsibility of the ultrasound system administrator at your facility to observe local security policies and implement security settings and policies for the ultrasound system.

You can lock the screen to prevent unauthorized access to patient information. An active screen saver display also locks the screen. A locked screen requires a user name and password to access the system.

See also: For information about the screen saver, refer to Chapter 1 in this manual.

The ultrasound system administrator also assigns permissions for each user account. Permissions establish data access rights (either no access or full access) that determine whether a user can view studies stored on the system. To understand the permissions assigned to your user account, or if you do not have a user name and case-sensitive password, contact your ultrasound system administrator.

Use the configuration settings to review the permissions for the currently logged-in user account.

Logging In to the Ultrasound System

All users are required to log in to the ultrasound system prior to use.

To log in to the ultrasound system:

1. Enter your user name and password on the login screen.
2. Click **OK**.

To	Do This
Prevent unauthorized access to the ultrasound system	1. Tap the login button.  2. Select Lock System.  The ultrasound system displays the login screen.
Begin a study without a password	(Available only when enabled in the configuration settings) Prerequisite: Obtain the emergency access user name from your ultrasound system administrator. 1. Enter the emergency access user name on the login screen. 2. Click OK. The ultrasound system enables limited access to patient data.

To	Do This
Change your password	(Available only when enabled for the currently logged-in user account in the configuration settings) 1. Tap the login button.  2. Select Change Password.  3. Enter the password information on the change password screen. 4. Click OK.
Log in as a different user	1. Tap the login button.  2. Select Log in Different User.  3. Enter your user name and password. 4. Click OK.
Log off the ultrasound system	1. Tap the login button.  2. Select Log Off <User Name> (where User Name is your assigned user name).  The ultrasound system displays a confirmation message to end the current session and log off the system. 3. Click Yes.

Registering a Patient and Beginning a Study

You can register a patient before you begin an exam or start the exam without registering a patient and then update the patient information at any time before ending the exam.

See also: For information about the fields in the registration form, refer to Appendix B in this manual.

Preregistering a Patient

When the ultrasound system is disconnected from the network, you can preregister a patient and begin a study. Once the study is complete, the system saves the preregistered information in the local database.

To preregister a patient and begin a study:

1. Press **PATIENT**.
2. Tap **New Patient**, if necessary.
3. Enter the required patient information.
4. Click **Preregister**.

The system stores the information on the preregistered page until the study is complete.

To begin a study from an existing preregistered page:

1. Click the **Worklist** tab.
2. Select **Preregister** from the **Server** list.
3. Click the required patient data.
4. Edit the study information, if necessary.
5. Tap **Start study**.

Registering a New Patient and Beginning a Study

To register a patient and begin a study:

1. Press **PATIENT**.
2. Tap **New Patient**, if necessary.
3. Enter the required patient information.
4. Tap **Start Study**.

The ultrasound system exits patient registration and activates real-time imaging.

5. To change the transducer and exam:

- a. Tap the transducer name.
- b. Tap the exam name.

The transducer and imaging settings activate for the selected exam.

Selecting an Existing Patient and Beginning a Study

You can select an existing patient from the local database on the ultrasound system.

When you begin a new study for an existing patient, the system populates the patient registration form with the preexisting patient data.

See also: For information about importing patient studies to the local database, refer to Chapter 5 in the System Reference.

To select an existing patient and begin a new study:

1. Press **PATIENT**.
2. Select the storage location for patient studies.
 - Tap **Browser**.
 - Click the **Local Database** tab.

See also: For information about searching for a patient study, refer to Appendix B in this manual.

3. Click the required patient data.
4. Edit the study information, if necessary, for example, the name of the exam.
5. Tap **New Study**.
The ultrasound system exits patient registration and activates real-time imaging.
6. To exit without beginning the patient study, tap **Cancel**.

Registering a Patient from the Worklist Server

When you select a patient from the worklist server, you can edit the patient registration information and begin the patient study.

Use the configuration settings to configure the worklist server and to assign a default exam and measurement package to a worklist procedure description.

To select a patient from the worklist and begin a new study:

Prerequisite: Connect the ultrasound system to a worklist server configured for modality performed procedure step (MPPS).

1. Press **PATIENT**.
2. Click **Worklist**.

The ultrasound system displays a list of patient studies on the worklist server.

See also: For information about searching for a patient study, refer to Appendix B in this manual.

3. Click the required patient.
4. Edit the study information, if necessary.
5. Tap **Start study**.
6. To exit without beginning the patient study, tap **Cancel**.

Customizing the Patient Registration Form

You can customize the display of patient information on the patient registration form.

To customize the patient registration form:

1. Press **PATIENT**.
2. Tap **New Patient**, if necessary.
3. Click the additional options button in the upper right of the patient registration form.

4. Select the information to display in the patient registration form.
See also: For information about information on the patient registration form, refer to Appendix B in this manual.
5. To include all information, click **Select All**.
6. To cancel the selections and retain the patient name, patient ID, transducer, and exam information, click **Deselect All**.

Editing Patient Registration Information

You can edit patient registration information at any time.

To edit patient information:

1. Press **PATIENT**.
2. Tap **Browser**, if necessary.
3. Select the required patient study.
4. Tap **Edit Patient**.
5. Enter the required information in the registration form.
6. Click **OK**.
7. To discard the changes, click **Cancel**.

Beginning a Study without Registering the Patient

You can begin a study without registering the patient.

To begin a study without registering the patient:

1. Begin the required scan.
2. Press the control assigned to the image or clip store function.
 The system saves the patient registration information as an image, inserts an asterisk (*) for the patient name, and uses the current date and time to assign a patient identification number.

Restarting an Existing Study

You can restart an existing study to continue imaging a patient, acquire images and clips, perform measurements and update calculated results, add annotations, and edit the patient registration information.

You can restart an existing study in the local database within 24 hours of starting the study.

To restart an existing study:

Prerequisite: The patient study is stored on the local database on the ultrasound system.

1. Press **PATIENT**.
2. Tap **Browser**.
3. Click the required patient study.
4. Tap **Restart**.
5. Edit the patient information, if necessary.
6. Tap **Restart**.

Selecting or Customizing the Exam Type

An exam type consists of a set of optimized imaging settings and a compatible transducer. The system includes the following factory-defined exams:

- **Abd** (Abdominal)
- **Abd-Difficult** (Abdominal Difficult)
- **Renal**
- **Bowel**
- **FAST** (Focused Assessment with Sonography in Trauma)
- **Lung**
- **OB** (Obstetrics)
- **Early OB** (Early Obstetrics)
- **Fetal Echo**
- **OB(Adv)** (Advanced Obstetrics)
- **GYN** (Gynecology)
- **Pelvic Floor**
- **Breast**
- **Testis**
- **Penile**
- **Thyroid**
- **Ped Hip** (Pediatric Hip)
- **Ped Abd** (Pediatric Abdomen)
- **Neo Head** (Neonatal Head)
- **Carotid**
- **Venous**
- **Arterial**
- **TCD** (Transcranial)
- **Pelvis**
- **Prostate**
- **Spine**
- **MSK** (Musculoskeletal)
- **Nerve**
- **Digital**
- **IO Abd** (Intraoperative Abdomen)
- **Cardiac**
- **Ped Echo** (Pediatric Echo)
- **Neo Echo** (Neonatal Echo)
- **TEE** (Transesophageal Echocardiography)

Customizing Exams

You can create, edit, and delete user-defined exams. You can also change the default exam type assigned to the selected transducer.

Use the configuration settings to assign a measurement package, a library of text annotations, and a library of body markers for the exam.

To create a user-defined exam:

- During imaging, adjust the imaging parameters and settings required for the user-defined exam.

To	Do This
Create a new exam	1. Tap Add Exam. 2. Enter a name for the exam in Exam Name. 3. Select a measurement package, a library of text annotations, and a library of body markers. 4. To assign the created exam as the default for the selected transducer, select the Set as Default for Selected Transducer check box. 5. Click Create. The system activates the created exam.
Create a new exam based on an existing exam	1. Tap the name of the transducer. 2. Tap Edit. 3. Select the required exam and then enter the exam name in Custom Exam. 4. Click Create. The system updates the exam list. 5. Click Close. The system activates the created exam.

To edit a user-defined exam:

1. Tap the required transducer tab and then tap **Edit**.
2. To rename the exam:
 - a. Select the required exam and then enter the exam name in **Custom Exam**.
 - b. Click **Rename**.
3. To change the default exam for the selected transducer, select the required exam and then click **Set Default**.
4. To edit the list of exams:
 - Click **Add Blank** to insert a blank space before the selected exam.
 - Tap the up arrow button to move the selected exam up one position.



- Tap the down arrow button to move the selected exam down one position.



The system updates the list of exam selections on the touch screen.

5. To restore the settings, click **Restore Settings**.
6. To delete the exam, select the required exam and then click **Delete**.
Note: You must exit a user-defined exam before you can delete the exam.

Activating and Deactivating Operating Modes

When you begin an exam, the system activates 2D imaging. You can activate multiple operating modes, for example, 2D with Color or Doppler. When you activate a mode, the system assigns priority to the selected mode and displays touch screen controls for the priority mode. You can assign priority to another mode at any time.

To activate or deactivate an operating mode:

1. Press the required control on the control panel. For example, press **C** to activate color imaging.
Note: If the system is configured to display the Doppler or M-mode cursor before activating the Doppler spectrum or M-mode sweep, adjust the position of the cursor and then press the required control again. Or, press **UPDATE**.
2. To deactivate all operating modes except 2D-mode, press **2D**.
3. To assign priority to another operating mode during combined mode imaging, tap the required mode on the touch screen. For example, tap **Color** to assign priority to color imaging during 2D imaging with color.

Workflow Protocol

A workflow protocol is a pre-defined checklist that guides you through a clinical workflow. The protocol consists of one or more views. A protocol view is an image or clip that can also include measurements, a body marker, and annotations. When the image or clip is saved and the measurements are performed, a view is complete.

Use the configuration settings to create and customize user-defined workflow protocols.

To activate and select a workflow protocol:

1. Tap **Protocols** and then select the required protocol view on the touch screen.
The system displays a blue border around the active view.
2. Press the control assigned to the print and store function, for example, **PRINT STORE 1**.
If the protocol view includes a measurement, the system selects the required measurement label and activates the default measurement method for the label.
3. Position and anchor all measurement markers required for the measurement, if necessary.

To temporarily modify a current workflow protocol:

To	Do This
Pause the protocol view	▪ Tap Pause Protocol. Or, press SELECT. Images saved during a paused protocol are not saved as protocol views.
Resume a paused protocol	▪ Tap Resume Protocol. Or, press SELECT.
Copy a protocol view	▪ Tap Copy View.
Change the name of a protocol view	1. Tap Rename View. 2. Enter a name and then click OK.
Delete a protocol view	▪ Tap Delete View. Note: Deleted measurement will not be deleted from the report.
Change the progression for the views	▪ Tap Progression. – Automatic activates the next view in the protocol when the current view is complete. – Manual waits for you to select the next view in the protocol. You can capture multiple images for a single view.

Optimizing Images

During real-time imaging, you can optimize the image using the control panel and on-screen controls.

Use the configuration settings to customize the imaging parameters assigned to on-screen controls.

Adjusting the Settings for General Imaging Parameters

The following imaging parameters apply to all operating modes.

To	Do This
Adjust the digital gain compensation	Slide the DGC controls on the control panel to the required positions.
Adjust the gain for the priority mode	Rotate the required mode control on the control panel. For example, to adjust color gain, rotate C.
Adjust the transmit frequency for the active transducer	Rotate [Frequency]. A higher frequency typically increases resolution; a lower frequency improves penetration. Note: Available settings depend on the transducer.
Adjust transmit power for the active transducer	Rotate [Tx Power].
Adjust the imaging depth	Rotate or push DEPTH / FOCUS as specified in the configuration settings. The system displays the depth in centimeters (cm) below the image.
Adjust the focal zones	- Press DEPTH / FOCUS to activate the multiple transmit focal zones and then push or rotate the control to reposition the transmit focal zones in the image as specified in the configuration settings. - To adjust the distance between multiple focal zone markers, rotate [Focus Width].
Magnify an image	(Not available for a rotated image) - Rotate ZOOM to adjust the magnification. The system activates the zoom function, a Z displays below the image, and the magnification level displays in the imaging parameters. - To pan the magnified image, roll the trackball. - To cancel magnification, press ZOOM.

To	Do This
Magnify a high-definition (HD) region of interest	(Not available for a rotated image) 1. Press ZOOM. 2. Adjust the size or position of the ROI. a. Press a trackball key, if necessary. ▪ Position indicates the trackball is assigned to repositioning the ROI. A solid line outlines the ROI. ▪ Size indicates the trackball is assigned to resizing the ROI. A dotted line outlines the ROI. b. Roll the trackball to adjust the size or position. 3. Press ZOOM again. When the image is magnified, the system displays HD Z next to the depth indicator on the image screen. 4. Adjust the size or position of the high-definition ROI, if necessary. 5. To cancel magnification, press ZOOM.
Pan a magnified image	▪ Roll the trackball to pan the image in the required direction.
Increase the image size to hide the left panel, thumbnail panel, and trackball legend	▪ Tap Extensive View.
Increase the image size to hide the left panel, thumbnail panel, trackball legend, and patient banner	▪ Tap Full Screen. Note: The captured images or clips in Full Screen do not include the patient banner.

2D Imaging

During 2D imaging with a curved or phased array transducer, the system displays an adjustable field of view (FOV).

Optimizing 2D Images

You can use an image preset to apply predefined image optimization settings for 2D. You can also adjust imaging settings individually to refine the optimization.

To	Do This
Resize or reposition of the field of view (FOV)	(Available only for phased and curved array transducers) 1. Press SET, if necessary. – Position indicates the trackball is assigned to repositioning the FOV. A dotted line outlines the FOV. – Size indicates the trackball is assigned to resizing the FOV. A dotted line outlines the FOV. 2. Roll the trackball to adjust the size or position and then press SET.
Reduce or enlarge the size of the 2D image	1. Tap the required button in FOV Size to adjust the size of the 2D image. – To reduce the image size, tap the left arrow.  – To enlarge the image size, tap the right arrow.  2. To restore the image size, tap Full FOV.

To	Do This
Extend the field of view to 224°	(Available only for the 9MC3 transducer) ▪ Tap Wide FOV.
Adjust the overall contrast resolution (dynamic range) of an image	▪ Rotate [DR].
Select a processing curve that assigns echo amplitudes to gray levels	▪ Rotate [Map].
Change the color tint applied to the 2D image	▪ Rotate [Tint].
Create a visible smoothing effect by maintaining lines of image data for each frame of image data	▪ Tap Persist and then tap the required setting. A higher setting includes more frames in the image data.
Distinguish the contours of a structure during real-time imaging using edge enhancement	▪ Tap Edge and then tap the required setting.
Activate Harmonics to enhance visualization, improve image contrast and spatial resolution, or reduce noise	▪ Tap Harmonics.
Adjust the balance between the image line density and the frame rate	▪ Rotate [Line Density].
Activate spatial compounding	▪ Rotate [Compounding].
Reduce speckle and enhance contrast using the UltraArt universal image processing feature	▪ Rotate [UltraArt] to select the required setting: – Off disables the function. – Low applies the low level of speckle reduction. – Med (Medium) applies the medium level of speckle reduction. – High applies the high level of speckle reduction.
Reverse the orientation of the image	▪ Tap L/R Flip to reverse the horizontal (left/right) orientation. ▪ Tap U/D Flip to reverse the vertical (up/down) orientation. Note: If an image is rotated in 90-degree increments, you cannot reverse the vertical orientation.
Rotate the image	(Available only for a single 2D-image acquired with a linear array transducer) (Not available for a magnified image or a frozen image with HD zoom) 1. Tap Rotate Image. 2. Tap the required setting: – 90° rotates the image clockwise in 90-degree increments. – -90° rotates the image counterclockwise in 90-degree increments. – 0° resets the image to the starting position. Note: A rotated image cannot be scaled to fit the image area.

To	Do This
Enable or disable the display of one real-time and one frozen image in a side-by-side format	(Also available when color is activated) 1. Press DUAL to display an image on the left side of the screen. 2. Press DUAL again to display an image on the right side of the screen. 3. To change the active image, press UPDATE. 4. To exit dual display format, tap Exit Dual. 5. To exit dual display format and activate the real-time imaging, press 2D.
Display two real-time images from the same acquisition on the image screen	Use the configuration settings to specify the dual display format. ▪ Tap Live Dual. Note: During color imaging, the color ROI is displayed only on the right side of the screen.
Display four separately acquired images	1. Tap Quad. 2. To continue to display images, repeatedly tap Quad.
Enable or disable a trapezoid image formats	(Available only for linear array transducers) ▪ Tap Trapezoid.
Adjust the angle of the field of view	(Available only for linear array transducers) ▪ Push STEER.
Optimize the uniformity of image brightness in the field of view	1. Press the tissue equalization control.  An icon on the image screen indicates the feature is active.  2. To update the image quality, rotate [TEQ Gain]. 3. To exit the feature, press the control again.
Decrease artifacts in a 2D-mode image	1. Tap Clarify. 2. To adjust the clarify level, tap Clarify Levels.
Display a center line on the image	(Available only for linear and curved array transducers) ▪ Tap Center Line. A vertical dotted line displays in the middle of the image.
Display grid lines on the image	(Available only for linear and curved array transducers) ▪ Tap Grid. Horizontal and vertical grid display as dotted lines. The distance between markers indicates the image depth.
Display the ruler	(Available only for linear and curved array transducers) ▪ Tap Ruler. A ruler displays with a triangle icon at the top center of the image. The distance between markers indicates the image depth.

To	Do This
Navigate the high-definition (HD) region of interest while referencing the full field of view	Prerequisite: Enable zoom navigation. Use the configuration settings to enable zoom navigation. 1. Press ZOOM twice to activate HD zoom. Zoom navigation displays a representation of the full field of view, an ROI at the anatomical point of magnification, and a transducer orientation indicator. 2. Adjust the size or position of the ROI using zoom navigation. a. Press a trackball key, if necessary. ▪ Position indicates the trackball is assigned to repositioning the ROI. A solid line outlines the ROI. ▪ Size indicates the trackball is assigned to resizing the ROI. A dotted line outlines the ROI. b. Roll the trackball to adjust the size or position. 3. To cancel magnification, press ZOOM.

Color Imaging

When you activate color imaging, the system displays an adjustable color region of interest (ROI). The system supports the following functions.

- Color imaging applies color to the velocities and directions of blood flow within the ROI.
- Power imaging applies color to the energy measurement generated by the blood flow within the ROI.
- Doppler Tissue Imaging (DTI) applies color to the velocities of tissue in the ROI.
- Doppler Tissue Energy (DTE) applies color to the Doppler signal energy returning from the tissue.

To	Do This
Resize or reposition the color region of interest (ROI)	1. Press SET, if necessary. – Position indicates the trackball is assigned to repositioning the color ROI. A solid line outlines the color ROI. – Size indicates the trackball is assigned to resizing the color ROI. A dotted line outlines the color ROI. 2. Roll the trackball to adjust the size or position.

Acoustic Streaming with Cysts

⚠ WARNING: Acoustic streaming may occur in fluid filled cysts under certain conditions, for example, when it is necessary to repetitively image a fluid filled cyst with long ultrasound pulses. To verify that flow inside of a cyst is actually acoustic streaming, use spectral Doppler with a large gate and low scale and filter settings and move the sample volume around the inside of the cyst.

Acoustic streaming is a well-known physical phenomenon that can occur while imaging fluid filled cysts. Under typical imaging conditions, the color flow and Doppler scale and filter settings are set too high to see the low velocity acoustic streaming, however, physical artifacts may be generated under the following conditions:

- Long ultrasound pulses are repetitively used to image a fluid filled cyst about 1 centimeter or larger.
- The fluid in the cyst has low viscosity.
- The center of the cyst is located near the elevation focus of the transducer.
- The PRF setting is low. (Acoustic streaming flow is typically visible between 1 to 2 centimeters per second.)
- The wall filter setting is low. (Acoustic streaming velocities are usually set by filters.)

You can confirm the flow inside of a cyst is acoustic streaming by using a large gate and low scale and filter settings and then moving the sample volume inside the cyst.

Note: Flow away from the transducer is visible at multiple locations throughout the center of the cyst if it is acoustic streaming. It is more difficult to generate acoustic streaming near the edge of the cyst.

Optimizing Color Images

To	Do This
Activate color imaging	▪ Tap Color .
Activate power imaging	▪ Tap Power .
Activate Doppler Tissue Imaging or Doppler Tissue Energy	▪ Tap DTI or DTE .
Extend the range of displayed velocities in one direction by repositioning the color baseline	▪ Rotate [Baseline] .
Adjust the range of displayed color velocities by adjusting the scale factor of pulse repetition frequency (PRF)	▪ Rotate [Scale] .
Adjust the angle of the color region of interest	(Available only for linear array transducers) ▪ Push STEER .


To	Do This
Select an image preset to optimize hemodynamic flow conditions	▪ Rotate [Flow] to select the required setting: – Low provides maximum sensitivity to low velocity flows by incorporating the lowest possible filter settings and lower pulse repetition frequency. You may experience increased flash (motion artifacts). – Med (Medium) produces an optimal balance between flash suppression and maximum sensitivity by using an adaptive wall filter. – High optimizes the system for the high arterial flow common to pulsatile vessels and stenotic conditions.
Reduce clutter by adjusting the level of motion discrimination during real-time imaging	▪ Rotate [Filter].
Adjust the time that the colors corresponding to blood flow velocity remain in the region of interest before decaying or being replaced by another color	▪ Rotate [Persist].
Adjust the threshold for overlaying color pixel information on the 2D image	▪ Rotate [Priority].
Display or hide color information in the region of interest	▪ Tap Color Display.
Adjusts the balance between the image line density (resolution) and the frame rate (speed)	▪ Rotate [Line Density]. A color bar indicates the range of colors for the selected map.
Assign color to velocity or energy information	▪ Rotate [Map].
Adjust the level of spatial (both axial and lateral) averaging used to smooth the flow pattern display	▪ Rotate [Smooth].
Invert the colors depicting forward and reverse flow in the color ROI and in the color bar	▪ Tap Invert.
Designate or tag a specific blood flow velocity or a range of blood flow velocities in real-time or frozen image or during CINE review	▪ Tap VelTag and then select the required settings.
Display energy generated by the direction of blood flow	(Available only during power imaging) ▪ Tap Dir Power.
Adjust the blending of color and 2D data on the image	▪ Rotate [α-Blend].

To	Do This
Navigate the high-definition (HD) region of interest while referencing the full field of view	Prerequisite: Enable zoom navigation. Use the configuration settings to enable zoom navigation. 1. Press ZOOM twice to activate HD zoom. Zoom navigation displays a representation of the full field of view, an ROI at the anatomical point of magnification, and a transducer orientation indicator. 2. Adjust the size or position of the ROI using zoom navigation. a. Press a trackball key, if necessary. ▪ Position indicates the trackball is assigned to repositioning the ROI. A solid line outlines the ROI. ▪ Size indicates the trackball is assigned to resizing the ROI. A dotted line outlines the ROI. b. Roll the trackball to adjust the size or position. 3. To cancel magnification, press ZOOM.

Pulsed Wave Doppler

When you activate pulsed wave Doppler, the system displays the Doppler cursor and gate on the 2D image before activating the Doppler spectrum. You can adjust the position of the cursor and the size of the gate before you activate the Doppler spectrum.

Use the configuration settings to define your pulsed wave Doppler preferences, for example, bypassing the cursor adjustment step before displaying the Doppler spectrum, selecting the units (Hz or cm/s) displayed in the Doppler spectrum, and selecting the display format for the Doppler spectrum and 2D image.

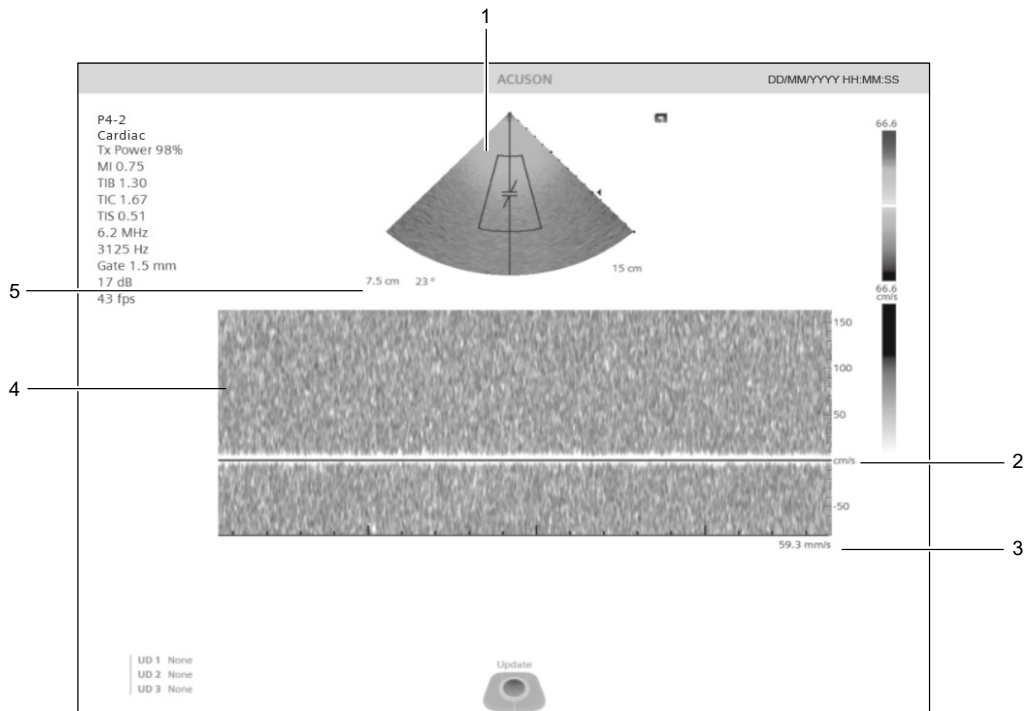

Example of image screen with pulsed wave Doppler activated.

- 1 2D image
- 2 Doppler units (frequency or velocity)
- 3 Sweep Speed
- 4 Doppler spectrum
- 5 Flow angle value in degrees and depth of the Doppler gate in centimeters

Optimizing the Doppler Spectrum

To	Do This
Activate or deactivate steerable continuous wave Doppler	(Available only for phased array transducers) Press CW.
Reposition the Doppler cursor	Roll the trackball.
Resize the Doppler gate	Rotate [Gate].
Pause the Doppler spectrum and refresh the 2D image	Press UPDATE.
Adjust the overall contrast resolution (dynamic range) of an image in 5 dB increments	Rotate [DR].
Display only the Doppler spectrum and hide the 2D image	Tap Full D.
Display both the 2D image and the Doppler spectrum	- Tap the split screen button.  - Tap the required format. 
Correct the flow angle	- To adjust the flow angle in 60° increments (60°, 0°, and -60°), press ANGLE. - To adjust the flow angle in 1° increments, rotate ANGLE.
Extend the range of displayed velocities by repositioning the Doppler baseline	Rotate [Baseline].
Adjust the scale factor of pulse repetition frequency (PRF) to improve flow sensitivity or to control aliasing	Rotate [Scale].
Adjust the interval for refreshing the 2D image	Tap Update Rate and then tap the required setting.
Adjust the scrolling speed of the Doppler spectrum	Rotate [Sweep].
Select a processing curve that assigns echo amplitudes to gray levels	Rotate [Map].
Change the color of the Doppler spectrum by adding blue, red, yellow, or green to the gray in a gray map	Rotate [Tint].
Reject signals below the selected frequency setting, for example, signals generally caused by tissue clutter	Rotate [Filter].
Enable or disable the simultaneous display of the 2D and color image and Doppler spectrum in real time	Tap Simult.

To	Do This
Reverse the Doppler spectrum relative to the baseline and invert the Doppler audio	▪ Tap Invert.
Adjust the Doppler audio volume	▪ Rotate [Volume].
Enable or disable the audio volume	▪ Tap Mute.
View Doppler statistics on a real-time spectrum	1. Tap Auto Trace. 2. Tap Trace Direction and then select the method for analyzing the Doppler spectrum. – Above displays statistics for Doppler data above the baseline. – Below displays statistics for Doppler data below the baseline. – Both displays statistics for Doppler data above and below the baseline. The system positions the peak systole at the maximum velocity of the Doppler spectrum and the end diastole at the end of the heart cycle. 3. To adjust the sensitivity level of the detection, tap Sensitivity. 4. To display the velocities or frequencies on the spectrum as a mean graphical trace, tap Mean Trace. The blue trace indicates a mean velocity or frequency. The yellow trace indicates the maximum value.
View Doppler statistics on a frozen spectrum or a spectrum during review	Prerequisite: Freeze the Doppler spectrum or view the Doppler spectrum during review and then activate the measurement function. 1. Tap PI Auto. 2. Position and anchor the beginning and ending markers on the spectrum at the peak systole and end diastole. 3. To change the method for analyzing the Doppler spectrum, tap Trace Direction and then select the required method: – Above displays statistics for Doppler data above the baseline. – Below displays statistics for Doppler data below the baseline. – Both displays statistics for Doppler data above and below the baseline. 4. To adjust the sensitivity level of the detection, rotate [Sensitivity]. 5. To display the velocities or frequencies on the spectrum as a mean graphical trace with the maximum graphical trace, press [Mean Trace]. The blue trace indicates a mean velocity or frequency. The yellow trace indicates the maximum value.
Adjust the angle of the Doppler gate left or right	(Available only for linear array transducers) ▪ Push STEER.
Activate Doppler Tissue Imaging	(Available only for phased array transducers) ▪ Tap DTI.

M-mode

When you activate M-mode, the system displays the M-mode cursor on the 2D image before activating the M-mode sweep. You can adjust the position of the cursor before you activate the M-mode sweep.

Use the configuration settings to define your M-mode preferences, for example, bypassing the cursor adjustment step before displaying the M-mode sweep and selecting the display format for the M-mode sweep and 2D image.

To position the M-mode cursor:

- Use the trackball to position the cursor in the required location.

Optimizing the M-mode Sweep

To	Do This
Adjust the overall contrast resolution (dynamic range) of an image in 5 dB increments	▪ Rotate [DR].
Adjust the angle of the M-mode cursor	(Available only for linear array transducers) ▪ Push STEER.
Display or hide the 2D image	▪ Tap Full M.
Display both the 2D image and the M-mode sweep	1. Tap the split screen button.  2. Tap the required format. 
Adjust the scrolling speed of the M-mode sweep	▪ Rotate [Sweep].
Distinguish the contours of a structure during real-time imaging using edge enhancement	▪ Rotate [Edge].
Select a processing curve that assigns echo amplitudes to gray levels	▪ Rotate [Map].
Change the color of the M-mode sweep by adding blue, red, yellow, or green to the gray in a gray map	▪ Rotate [Tint].

Examples of Imaging Settings

The image screen displays exam and imaging settings for the priority operating mode.

Settings for All Modes

Example	Description
5C1	Transducer
Abdomen	Exam type
Tx Power 100%	Transmit power, a percentage of the maximum acoustic intensity and mechanical index for the active transducer
MI 0.78	Mechanical index
TIS 1.7	Soft tissue thermal index
TIC 1.7	Cranial bone thermal index
TIB 1.7	Bone tissue thermal index

2D Settings

Example	Description
Compound	Compounding imaging
UA 1	UltraArt universal image processing
DR 75 dB	Dynamic range
3.3 MHz	Frequency Note: An H displays prior to the frequency during Harmonics imaging.
29 dB	2D gain
36 fps	Frames per second
Zoom 1.1	Magnification level

Color Settings

Example	Description
2.3 MHz	Frequency
1220 Hz	Pulse Repetition Frequency (PRF)
Flow Med	Color flow image preset
5 dB	Color gain
8 fps	Frames per second

Pulsed Wave Doppler Settings

Example	Description
2.3 MHz	Frequency
1775 Hz	Pulse Repetition Frequency (PRF)
Gate 1.5 mm	Size of the Doppler gate
11 dB	Doppler gain
8 fps	Frames per second

M-mode Settings

Example	Description
2.3 MHz	Frequency
30 dB	M-mode gain
24 fps	Frames per second

General Contrast Agent Imaging Settings

Example	Description
Contrast	Imaging feature
1.8 MHz	Frequency
7 dB	Gain during contrast agent imaging
DR 60 dB	Dynamic range
P 2	Persistence, a visible smoothing effect
Peak Hold	Peak hold, the maximum intensity projection


CINE

CINE continuously stores recently acquired data in a memory buffer in all operating modes. You can view data stored in CINE memory in continuous playback or frame-by-frame review. You can also scroll through the images and add measurements and annotations to any frame.

CINE memory is fixed. However, the length of playback varies based on the complexity of the images. The simpler the image, the more images can be retained.

During CINE playback, the system displays a CINE bar representing the status of the CINE memory buffer. In dual display format, the system displays a CINE bar below each image.

Use the configuration settings to assign CINE playback to a soft key.

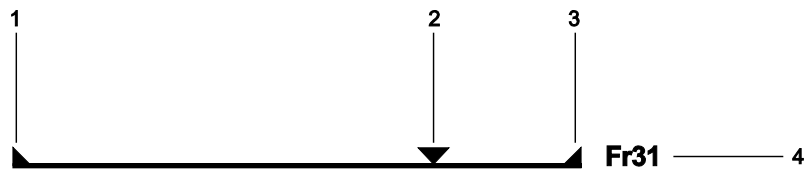

Example of a CINE bar.

- 1 Indicates the beginning of the CINE data. Reposition this marker to edit the length of the CINE data.
- 2 Indicates CINE playback. The marker depicts the location of the current frame within the loop of CINE data.
- 3 Indicates the end of the CINE data. Reposition this marker to edit the length of the CINE data.
- 4 Indicates the number of frames played back.

To review CINE data:

1. Press **FREEZE**.
The system displays the CINE bar.
2. To start or stop CINE playback, tap the play or pause button.
3. To adjust the speed of CINE playback, rotate **[Cine Speed]**.
4. To adjust the length of the clip, roll the trackball and then tap **Set Start Marker** or **Set End Marker** to position the current frame marker at the required location on the image screen.
5. To restore the original position of the markers, tap **Reset Markers**.
6. (Available only during pulsed wave Doppler) To toggle the 2D image with the Doppler spectrum during playback, press **UPDATE**.
7. To exit CINE, press **FREEZE**.

Capturing Images and Clips

When you capture an image or clip, the system saves the data to the hard drive and displays a thumbnail of the image or clip in the thumbnail panel. You can use the review function to view saved images and clips, perform measurements, add annotations and body markers, and add images to reports.

To capture an image or clip:

- Press the control assigned to the image or clip store function (for example, **PRINT STORE 1**, **UD 1**, **CLIP STORE**, or the optional footswitch).

Reviewing Images and Clips

Use the review function to view, evaluate, measure, and annotate the results of the current examination or a previous examination saved to the hard disk.

Printing Images

Note: Refer to the manufacturer's user manual for the proper handling and operating instructions for the device connected to your system.

You can print images from a current examination, a patient report, CINE data, or an image retrieved from a previous study.

To print an image:

- Press the control assigned to the print function (for example, **PRINT STORE 1**).

Video Recording a Study

Note: Refer to the manufacturer's user manual for the proper handling and operating instructions for the device connected to your system.

The ultrasound system does not have controls for video playback.

To record a study:

Prerequisite: Use the configuration settings to assign the video recording function to the user-defined controls on the control panel.

- Use the controls on the video recording device. Or, press the control assigned to the video recording function (for example, **UD 1** or **UD 2**) to record and pause video recording.

To stop recording:

- Use the controls on the video recording device. Or, end the current patient examination.

When recording on a video recording device is in progress or paused, the ultrasound system displays a video icon in the system status bar.

Icon	Description
	Indicates data recording is in progress.
	Indicates data recording is paused.

Completing an Examination

You can end an examination from the review screen or the patient registration screen.

To end the current patient examination, choose a method:

- Tap **End Exam**.

5 Examination Review


Reviewing Current and Previous Exams	2
Activating Review	2
Playing and Editing Clips.....	4
Exporting Images and Clips.....	5

Reviewing Current and Previous Exams

Use the review function to view and evaluate the results of the current examination or a previous examination saved in the local database (system hard disk).

Activating Review

When you activate review, the system displays all images and clips saved with the study. You can also view archived reports during review.

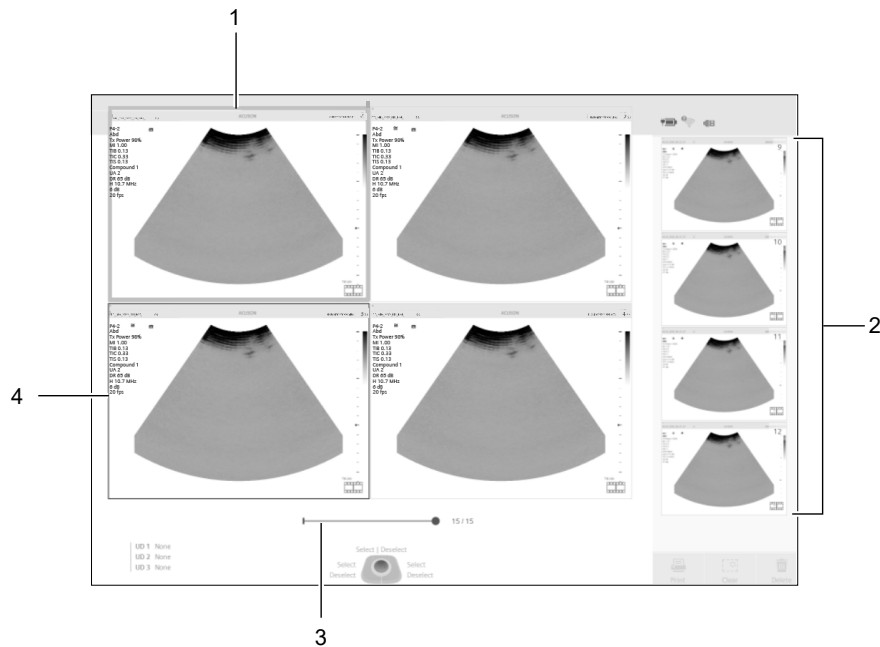

Example of review screen.

- 1 Active image.
The system displays a wide transparent white border around the currently active image.
- 2 Thumbnail panel
- 3 CINE bar
- 4 Selected image.
The system displays a white border around a selected image.

To activate review:

1. To activate review during an exam, press **REVIEW**. Or, double-click a thumbnail.
2. To activate review of a saved exam:
 - a. Tap **Browser**.
 - b. Select a patient from the patient list and then select a study from the study list.
 - c. Double-click the selected study.
3. To change the layout of the images:

To	Do This
Change the display format of the saved images and clips	Selects a display format on the touch screen: ▪ 1x1 ▪ 1x2 ▪ 2x2 ▪ 3x3 ▪ 4x4 ▪ 5x5
Scroll through multiple pages of images or clips	▪ Rotate [Page].
View the active image or thumbnail in a full-screen layout	▪ Double-click the active image or required thumbnail.
Select an image on the review screen or thumbnail panel	1. Click the image or required thumbnail. 2. To cancel the selection, click the image again. 3. To select all images on the review screen, tap Select all. 4. To cancel the selection of all images, tap Deselect All.
Print selected images on the review screen or thumbnail panel	1. Click Print. 2. Select the required printer and then click OK.
Delete selected images on the review screen or thumbnail panel	1. Tap Delete. 2. Tap OK to confirm the deletion. 3. To continue without deleting the images, click Cancel.
View the selected images on the review screen	1. Tap Show Selected. 2. To cancel the selection of selected images, tap Show All.
View the slide show for the selected images	1. Tap Slide Show. 2. To stop the slide show, click Stop. 3. To continue playing the slide show, click Start. 4. To exit the slide show and activate review, click Exit.
View a report from a saved exam	1. Tap Archived Report. 2. To exit the report, tap Close.
Activate the arterial health package (AHP)	1. Tap AHP. 2. To exit the screen, tap Exit.
View the carotid intima-media thickness (CIMT) report	1. Tap AHP. 2. Click CIMT Report.
Exit review and activate real-time imaging	▪ Press REVIEW. Or, press 2D.
Exit review and end the exam	▪ Tap End Exam.

Playing and Editing Clips

When the active image is a clip, the system activates the clip playback controls.

Use the configuration settings to configure your preferences for clip playback.

To	Do This
Play or pause playback of one or more clips	- To start playback, tap the play button.  - To pause playback, tap the pause button. 
Review the first or end CINE frame	- To review the first CINE frame, tap the first frame button.  - To review the end CINE frame, tap the end frame button. 
View a clip frame-by-frame	Rotate [Cine Frame] to display the next or previous frame in the clip.
Adjust the playback speed	Rotate [Cine Speed] to select the required speed. Note: Playback speed is a multiple of the speed at acquisition. For example, if the clips were acquired at 30 fps (frames per second), specify a playback speed of 1/2 to play back the clips at 15 fps.
Change the current frame	1. Pause clip playback. 2. Drag the current frame marker to the new location. Or, rotate [Cine Frame].
Change the length of a clip	1. To change the starting frame of the clip, rotate [Cine Start Marker]. 2. To change the ending frame of the clip, rotate [Cine End Marker]. 3. To restore the original clip length, tap Reset Marker.
Capture a clip	Tap Save.
Store the current frame of a clip	1. Pause playback of the clip. 2. Press the control on the control panel assigned to the image or clip store function.

Exporting Images and Clips

You can export images and clips to external media, for example, a DVD or USB storage device.

Note: You can also export images, clips, and complete studies from the patient browser.

To export images:

1. Select one or more images on the review screen and then tap **Export**.
The system displays a dialog box to select the format and destination for the exported images.
2. Select the required **Format**:
 - **jpeg / avi**
 - **DICOM**
3. To remove any patient-identifying information from the images:
 - Select the **De-identify** check box.

Note: You can enter text to replace **Anonymous** as the patient name.
4. Select the target location, for example, a USB port, for the exported images in the **Destination** list.
5. Enter the destination path for the export location in the **Directory** box.
6. Click **Export**.
The system displays progress information until the export is complete.
7. When the export process is complete, click **Close**.

6 Transducer Care and Accessories

Caring for Transducers	2
Methods for Cleaning and Disinfecting Transducers	4
Examples of Transducer Components	5
Immersion Levels	7
Cleaning, Disinfecting, and Sterilizing Transducers.....	8
Cleaning a Transducer at the Point of Use.....	9
Cleaning a Transducer	10
Disinfecting a Transducer.....	11
Disinfecting a Continuous Wave Transducer.....	12
Sterilizing a Transducer.....	12
List of Approved Cleaners	13
List of Approved Disinfectants	14
List of Approved Sterilization Systems and Sterilants	15
Transducer Accessories	16
Transducer Sheaths (Transducer Covers).....	17
Disposal.....	18
Storage.....	19
Brackets and Needle Guides	19

Caring for Transducers

See also: For information about cleaning and disinfecting transesophageal transducers, refer to Chapter D3 in the Clinical Applications Guide.

-  **WARNING:** To minimize the risk of cross-contamination and infectious diseases, clean and high-level disinfect endocavity transducers after each use; clean and sterilize intraoperative transducers after each use.
-  **WARNING:** To avoid electrical shock, you must visually inspect a transducer prior to use. Do not use a transducer with cracked, punctured, or discolored casing or a frayed cable.

Discoloration Exception: The use of approved cleaners and disinfectants may cause discoloration of the transducer. You can continue to use a transducer discolored due to the use of approved cleaners and disinfectants only.
-  **WARNING:** Prior to each use, inspect the endocavity or intraoperative transducer for signs of mechanical damage such as cracks, cuts, tears, perforations, or protrusions. Do not use the transducer if the transducer appears damaged in any way. Any damage could cut the patient and compromise the electrical safety of the transducer, causing possible patient or user injury. Contact your local Siemens Healthineers representative.
-  **WARNING:** During neurosurgical procedures, if a transducer becomes contaminated with tissue or fluids of a patient known to have Creutzfeldt-Jakob disease, the transducer should be destroyed, as it cannot be sterilized.
-  **WARNING:** When using an endocavity or intraoperative transducer with a CF type applied part, the patient leakage currents may be additive.
-  **Caution:** Transducers are sensitive instruments – irreparable damage may occur if they are dropped, knocked against other objects, cut, or punctured. Do not attempt to repair or alter any part of a transducer.
-  **Caution:** To avoid cable damage, do not roll the ultrasound system over transducer cables, and do not wrap the transducer cable around the transducer handle or the transducer connector.
-  **Caution:** To avoid damage to the transducer, do not use transducer sheaths containing an oil-based coating or petroleum- or mineral oil-based ultrasound coupling agents. Use only a water-based ultrasound coupling agent.
-  **Caution:** Follow all instructions provided by manufacturers of sterile goods (transducer sheaths) to ensure proper handling, storage, and cycling of all sterile goods.
-  **Caution:** To avoid damage caused by electrostatic discharge, do not touch the electrical component on the transducer connector.

Take extreme care when handling or storing transducers. They must not be dropped, jarred, or knocked against other objects. Do not allow transducers to come into contact with any sharp-edged or pointed object.

Protective Case

Due to the mechanical sensitivity of transducers, Siemens Healthineers recommends that you always use the transducer case when you ship a transducer or transport it from one place of examination to another. The case is specially designed to protect the sensitive parts of the transducer. Be sure that all parts of the transducer are properly placed inside the case before you close the lid.

Storage

Store transducers in a clean, dry, protected environment. Extreme temperatures or humidity may damage a transducer.

See also: For information about storage temperatures, refer to Appendix A in this manual.

Transducer Element Check

A transducer element check is available for evaluating the integrity of each element on a transducer. If you suspect transducer degradation due to damage, age, or another issue, contact your customer service representative. Based on the results of the transducer element check, your customer service representative will advise you of necessary actions.

Repair

Siemens Healthineers can only guarantee the safety and performance of transducers serviced by an authorized third-party. Contact your local Siemens Healthineers service representative when your transducer appears to be damaged or malfunctions in any way.

Methods for Cleaning and Disinfecting Transducers

You must observe the required methods for cleaning and disinfecting or cleaning and sterilization before the first use and after each use of a transducer.

Siemens Healthineers defines the methods for cleaning and disinfecting or cleaning and sterilization by classifying the use of the transducer according to the Spaulding Classification scheme observed by the United States Food and Drug Administration (FDA).

Example of Transducers	Use	Minimum Required Method	Classification
▪ Intraoperative ▪ Laparoscopic	Normally sterile tissue or body space	Cleaning and sterilization	Critical
▪ Endocavity ▪ Transesophageal	Intact mucous membranes	Cleaning and high-level disinfection	Semi-critical
▪ Non-invasive	Non-intact skin	Cleaning and high-level disinfection	Semi-critical
▪ Non-invasive	Intact skin	Cleaning and low-level disinfection	Non-critical

Important Information

The ultrasound system does not support transducers for the following uses:

- ACUSON Juniper system
 - Laparoscopic
- ACUSON Juniper Select system
 - Laparoscopic
 - Intraoperative
 - Transesophageal


Reprocessing Methods

The transducer cleaning and disinfecting instructions in this chapter were validated for efficacy with the following reprocessing methods:


- Cleaning reprocessing methods:
 - Spraying or wiping a transducer
 - Immersing a transducer
- Disinfecting reprocessing methods:
 - Spraying or wiping a transducer
 - Immersing a transducer
 - Use of an automated reprocessing device
- Sterilization reprocessing methods:
 - Use of an automated sterilization system

See also: Contact your customer service representative if you require additional information about reprocessing a transducer.


Examples of Transducer Components

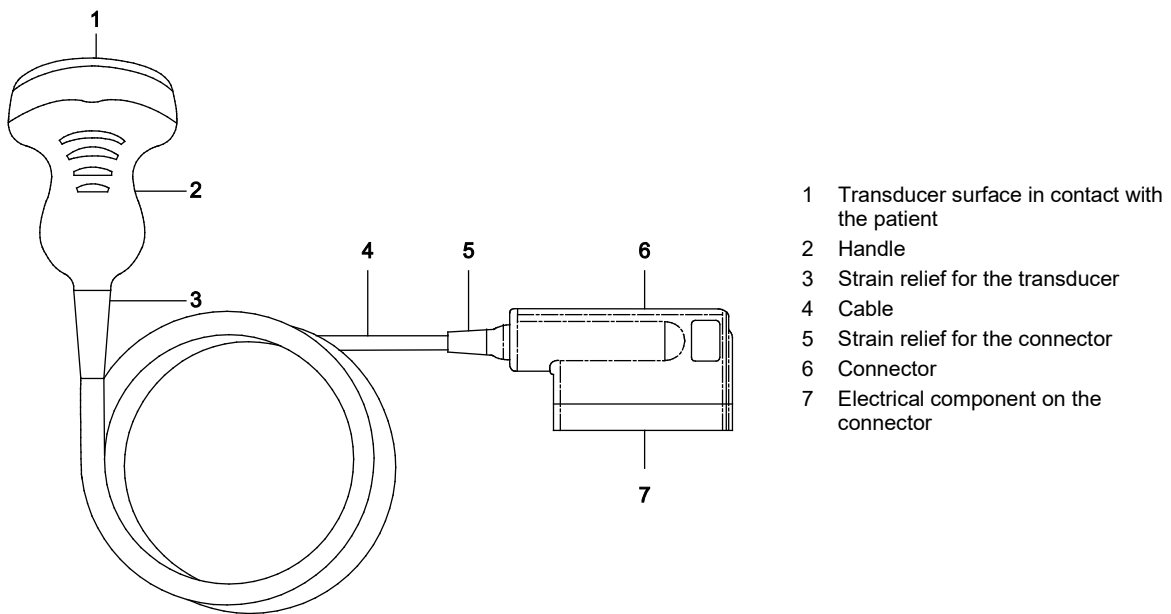

Example of transducer components.

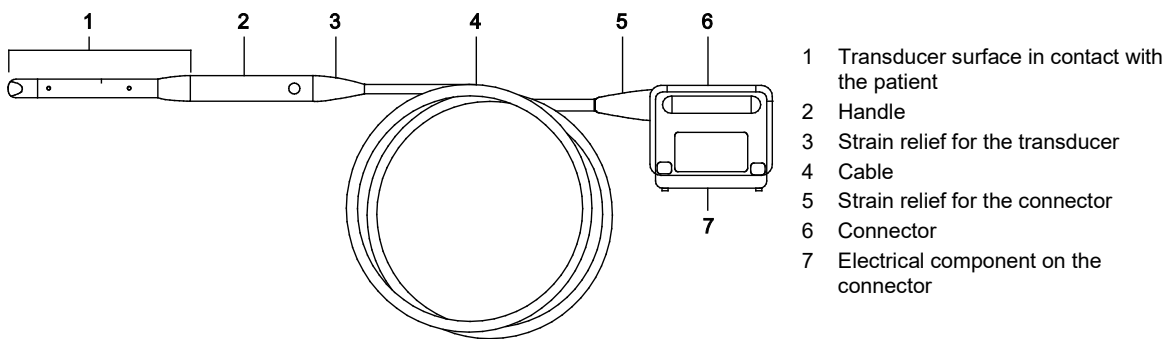

Example of endocavity transducer components.

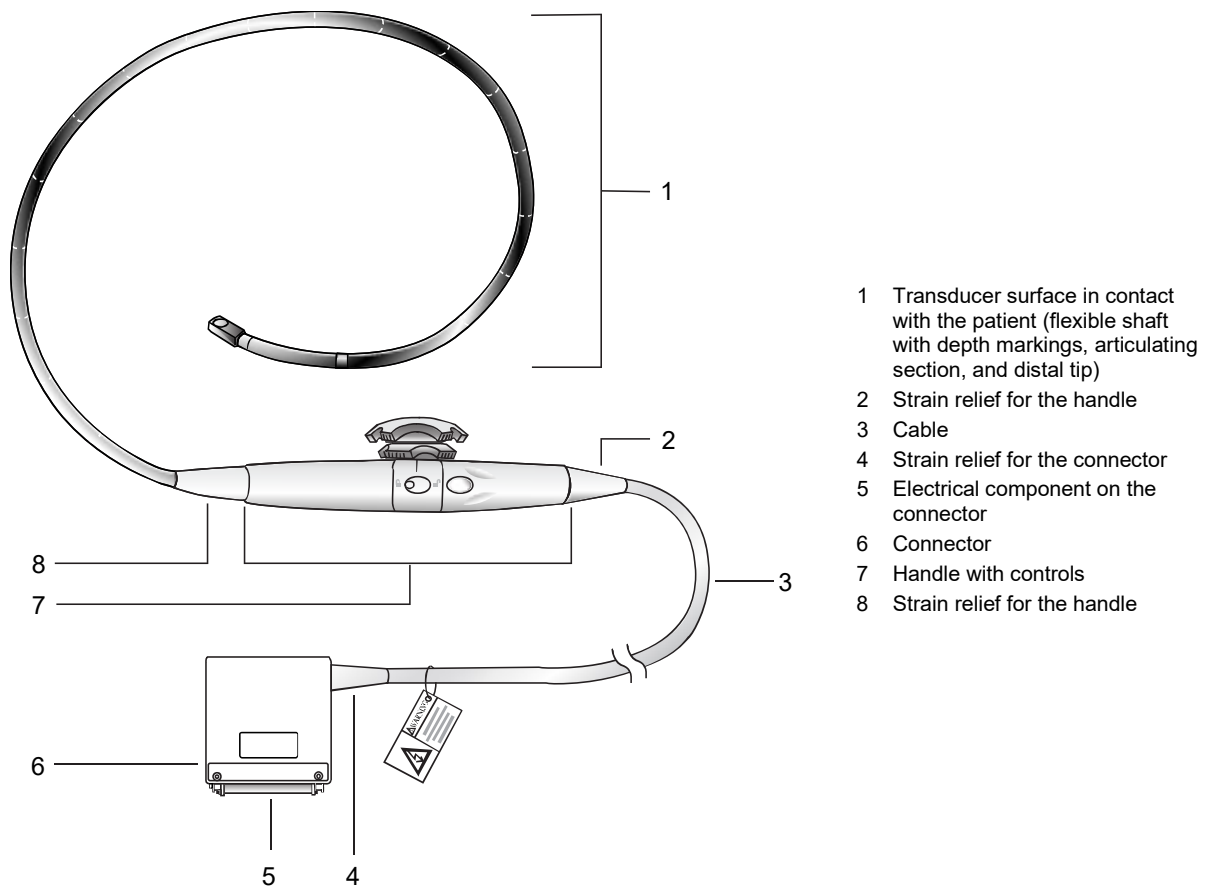


Example of transesophageal transducer components.

Immersion Levels

⚠ Caution: To avoid damage to the transducer, observe the immersion levels indicated for each transducer type. Transducers with the protection level IPX8 are indicated by the presence of the "IPX8" symbol on the connector of the transducer.

⚠ Caution: Do not immerse the label located on the cable of the CW transducer.

Note: Transducers meet Ingress Protection level IPX8 of EN 60529 and IEC 60529 to the depth of the immersion line shown in the illustration.

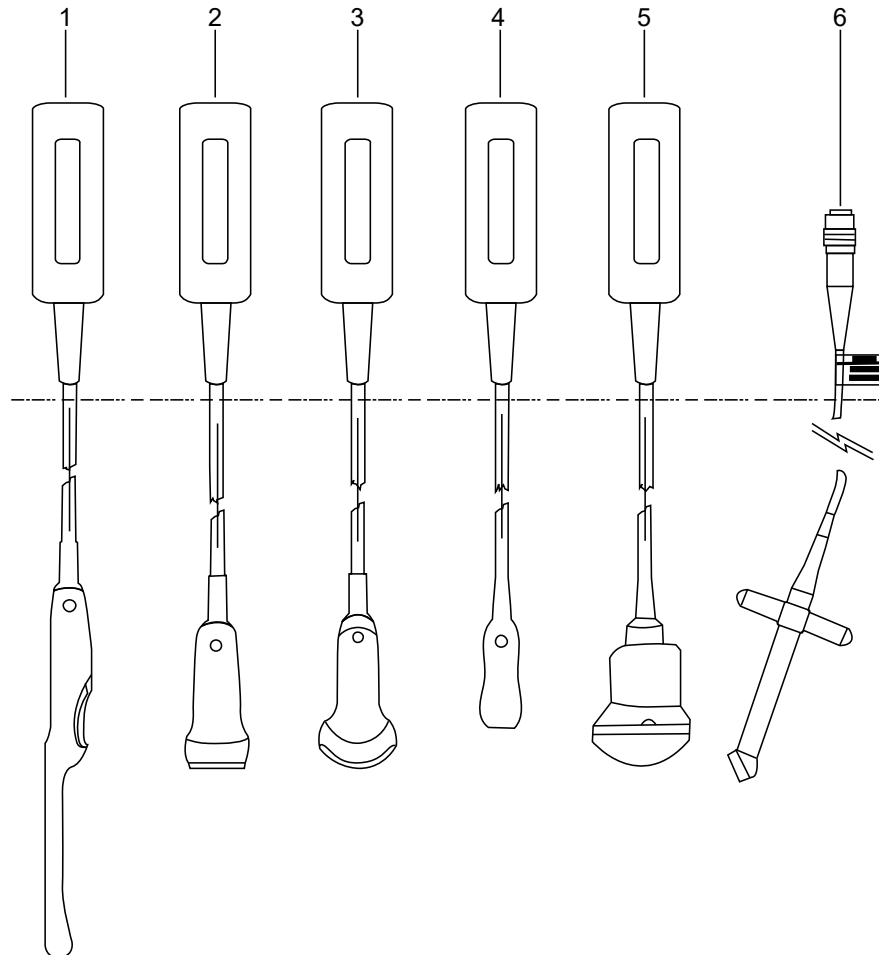

Example of immersion levels.

- 1 Endocavity
- 2 Linear
- 3 Curved
- 4 Phased
- 5 *fourSight* 4D
- 6 Continuous Wave (CW)

See also: For transesophageal transducer immersion levels, refer to Chapter D3 in the Clinical Applications Guide.

Cleaning, Disinfecting, and Sterilizing Transducers

-  **WARNING:** To avoid electrical shock and damage to the system, disconnect the transducer prior to cleaning or disinfecting.
-  **WARNING:** Disinfectants and cleaning methods listed are recommended by Siemens Healthineers for compatibility with product materials, not for biological effectiveness. Refer to disinfectant label instructions for guidance on disinfection efficacy and appropriate clinical uses.
-  **Caution:** Do not sterilize transducers using hot steam, cold gas, or Ethylene Oxide (EO) methods. Before applying any other methods which might be recommended by manufacturers of sterilization equipment, please contact your Siemens Healthineers representative.
-  **Caution:** To avoid damage to the transducer, observe the immersion levels indicated for each transducer type. Do not immerse or allow the connector of a transducer or the strain relief on the connector to become wet.
-  **Caution:** The transducers have been designed and tested to be able to withstand disinfection as recommended by the manufacturer of the disinfectant product. Carefully follow the disinfectant manufacturer's instructions.
-  **Caution:** Do not use abrasive cleaners, organic solvents such as benzene, isopropyl alcohol, or phenol-based substances, cleaners, or disinfectants containing organic solvents to clean or disinfect transducers. These substances can damage the transducers.

Isopropyl Alcohol Exception: Immersing a transducer in isopropyl alcohol can damage the transducer. Only use approved cleaners and disinfectants on the transducer.

You can only use a clean, lint-free gauze pad moistened with isopropyl alcohol to clean the electrical component on the transducer connector.
-  **Caution:** Do not use an abrasive sponge or brush. These materials can damage the transducer.
-  **Caution:** The use of a spray cleaner or disinfectant may force fluid inside the transducer and the ultrasound system. To avoid damage, first remove the transducer from the system and then carefully spray the transducer. Never spray the transducer connector.
-  **Caution:** Prior to cleaning or disinfecting a transducer, you must disconnect and remove the transducer from the ultrasound system. Do not clean or disinfect the transducer on or near the ultrasound system. Unintended contact with approved transducer cleaners and disinfectants can result in severe damage to the ultrasound system components.

Required Task	Location in This Manual
Read and understand the biohazard considerations.	Chapter 2
Perform the steps in the daily checklist.	Chapter 2
Read and understand the cleaning and disinfecting information for the ultrasound system.	Chapter 2
Read and understand all transducer-specific warnings and cautions.	page 6-2
Disconnect the transducer.	Chapter 3
Observe the required method of cleaning and disinfecting for the transducer.	page 6-4
Wipe visible soil and contaminants from the transducer immediately after each use.*	page 6-9
Clean the transducer.*	page 6-10
Disinfect the transducer.*	page 6-11
Sterilize the transducer, if required.	page 6-12
Store the transducer.	page 6-3

*For information about cleaning, disinfecting, and storing transesophageal transducers, refer to Chapter D3 in the Clinical Applications Guide.

Cleaning a Transducer at the Point of Use

Cleaning a transducer immediately after each use prevents soil and contaminants from drying on the surface of the transducer and facilitates subsequent cleaning.

To immediately clean a transducer after each use:

1. Disconnect and remove the transducer from the ultrasound system.
2. Remove any transducer accessories from the transducer.
3. Immediately clean the transducer in close proximity to the point of use.
 - a. Moisten a clean gauze pad or lint-free cloth with water or an approved cleaner. Or, use an approved cleaning wipe.
See also: For a list of approved cleaners, refer to page 6-13.
 - b. Carefully wipe the transducer surface intended for contact with the patient. Avoid touching the electrical components on the connector.
4. Perform subsequent cleaning of the transducer according to the next procedure.

Cleaning a Transducer

Cleaning a transducer removes soil and contaminants from the surface of the transducer and prepares the transducer for disinfection or sterilization.

To clean a transducer:

Prerequisite: The transducer is disconnected and removed from the ultrasound system. Transducer accessories are removed from the transducer, including the transducer sheath.

1. Select an approved cleaner. Follow the manufacturer's instructions for concentration and temperature of the solution and contact time when cleaning the surface of the transducer.
See also: For a list of approved cleaners, refer to page 6-13.
2. Clean the transducer from the connector to the surface intended for contact with the patient using one of the following cleaning methods. Avoid touching the electrical components on the connector.
 - If immersion is required, keep the transducer connector and the connector strain relief dry while immersing the transducer in an approved cleaner to the level indicated on page 6-7.
 - If wiping is required, carefully wipe the transducer. Continue wiping the surfaces of the transducer until there is no visible residue.
 - If spraying is required, do not spray the transducer near the ultrasound system. Never spray the transducer connector.
3. Rinse the transducer according to the cleaner manufacturer's instructions. Keep the electrical components on the connector dry.
Note: If the cleaner manufacturer does not provide rinsing instructions, thoroughly rinse the transducer with deionized water at room temperature.
4. Disinfect the transducer or sterilize the transducer, as required.

To clean the electrical component on the transducer connector:

1. To remove any gel from the electrical component, carefully wipe the electrical component using a lint-free gauze pad moistened with water.
2. To remove any skin oils from the electrical component, carefully wipe the electrical component using a lint-free gauze pad moistened with isopropyl alcohol.
3. Use compressed air to remove particles from the electrical component.

Disinfecting a Transducer

Immediately after cleaning, disinfect a transducer to destroy pathogens and other microorganisms and to prepare the transducer for storage or the next clinical use.

To disinfect a transducer:

Prerequisite: Clean the transducer. Identify the required level of disinfection for the transducer.

1. Select an approved disinfectant for the required level of disinfection. Follow the disinfectant manufacturer's instructions for concentration and temperature of the solution and contact time when disinfecting the surface of the transducer.

See also: For a list of approved disinfectants, refer to page 6-14.

2. Disinfect the transducer from the connector to the surface intended for contact with the patient using one of the following disinfection methods. Avoid touching the electrical components on the connector.

See also: For information about disinfecting a continuous wave (CW) transducer with a T-bar, refer to page 6-12.

- If immersion is required, keep the transducer connector and connector strain relief dry while immersing the transducer in an approved disinfectant to the immersion level indicated on page 6-7.
- If wiping is required, carefully wipe the transducer.
- If spraying is required, carefully spray the transducer. Do not spray the transducer near the ultrasound system. Never spray the transducer connector.
- If a reprocessing device is required, use a device compliant with ISO 15883 and follow the reprocessing device manufacturer's instructions.

3. Rinse the transducer according to the disinfectant manufacturer's instructions. Keep the electrical components on the connector dry.

Note: If the disinfectant manufacturer does not provide rinsing instructions, thoroughly rinse the transducer with deionized water at room temperature.

4. Follow the disinfectant manufacturer's instructions for drying the transducer.

Note: If the disinfectant manufacturer does not provide drying instructions, dry the transducer with a clean, lint-free, soft cloth.

5. Store the transducer in a clean, dry, protected environment.

Disinfecting a Continuous Wave Transducer

A continuous wave (CW) transducer with a T-bar requires additional tasks for disinfection.

To prepare a CW transducer for disinfection:

1. Rotate the screw on the side of the T-bar counterclockwise and then set the screw aside.
2. Loosen and then slide the T-bar onto the transducer cable.
3. Perform the additional tasks when following the disinfection instructions on page 6-11.
 - Carefully wipe the T-bar and screw using a disinfectant wipe.
 - Use a circular motion when wiping the T-bar and screw.
 - Slide the wipe partially through the hole on the T-bar and then rotate the T-bar to disinfect the interior of the hole.
 - Slide a portion of the wipe into the threaded screw hole on the T-bar and rotate the wipe to clean the interior of the hole.
 - Carefully wipe the transducer, including the cable, using a second disinfectant wipe.
 - Allow all surfaces to remain wet with disinfectant for a minimum of two minutes. Reapply the disinfectant to the transducer, T-bar, and screw, if necessary.

Note: Follow the disinfectant manufacturer's instructions.

 - Rinse and then dry the transducer, T-bar, and screw.
4. Reassemble the transducer after all components are dry.
 - a. Carefully slide the T-bar into place on the transducer.
 - b. Rotate the screw clockwise into the T-bar to secure the T-bar to the transducer.

Sterilizing a Transducer

Sterilization is required for transducers in contact with normally sterile tissue or body space, for example, transducers used in an intraoperative or laparoscopic environment. Sterilize a transducer to destroy pathogens and other microorganisms and to prepare the transducer for the next clinical use.

To sterilize a transducer:

Prerequisite: Clean the transducer. Disinfect the transducer, if necessary.

1. Select an approved sterilization system and sterilant for the transducer.

See also: For a list of approved sterilization systems and sterilants, refer to page 6-15.
2. Ensure the transducer is clean and dry before placing into the sterilization system.
3. Carefully follow the sterilization system manufacturer's instructions.
4. Follow the sterilization system manufacturer's instructions for removing the transducer from the sterilization system.
 - Use sterile technique when removing the transducer from the sterilization system.
 - Examine the transducer for any signs of damage and sterilant residue.
5. Follow the sterilization system manufacturer's instructions for drying the transducer.

List of Approved Cleaners

See also: For an additional list of approved cleaners and disinfectants, refer to the Transducer Cleaners and Disinfectants Addendum. You can also access the most recent compatible cleaners and disinfectants on the Siemens Healthineers web site.

Note: An approved cleaner may discolor the transducer. There is no associated degradation of imaging performance or transducer reliability.

	Cline Universal Wipes & Spray	ENZOL	Prolystica 2X Concentrate Enzymatic Presoak and Cleaner	Protex ULTRA Disinfectant Wipes	Transeptic Cleaning Solution
5C1	✓	✓	✓	✓	✓
6C1	✓	✓	✓	✓	✓
7C2	✓	✓	✓	✓	✓
12L3	✓	✓	✓	✓	NA
11L4	✓	✓	✓	✓	✓
13L4	✓	✓	✓	✓	✓
14L4	✓	✓	✓	✓	✓
16L4	✓	✓	✓	✓	✓
18H5	NA	✓	✓	✓	NA
11M3	NA	✓	✓	✓	✓
L10-5v	NA	✓	✓	✓	✓
5P1	✓	✓	✓	✓	✓
8V4	✓	✓	✓	✓	✓
10V4	✓	✓	✓	✓	✓
P4-2	✓	✓	✓	✓	✓
9MC3	✓	✓	✓	✓	✓
10MC3	✓	✓	✓	✓	✓
10BCC3	NA	✓	✓	✓	✓
9VC2	NA	✓	✓	NA	✓
8VC3	NA	NA	NA	✓	NA
9VE3	NA	NA	NA	NA	✓
9VE4	NA	✓	✓	NA	✓
CW2	✓	✓	✓	✓	✓
CW5	✓	✓	✓	✓	✓
CW8	NA	✓	✓	NA	✓

✓ = approved

NA = not approved

List of Approved Disinfectants

The Federal Republic of Germany prohibits the use of wipes for disinfection of transesophageal transducers and any transducer used in a semi-critical procedure.

See also: For an additional list of approved cleaners and disinfectants, refer to the Transducer Cleaners and Disinfectants Addendum. You can also access the most recent compatible cleaners and disinfectants on the Siemens Healthineers web site.

Note: An approved disinfectant may discolor the transducer. There is no associated degradation of imaging performance or transducer reliability.

	CIDEX OPA	gigasept FF (new)	Revital-Ox RESERT HLD	Sani-Cloth AF3 Germicidal Disposable Wipes	Super Sani-Cloth Germicidal Disposable Wipe	trophon EPR, NanoNebulant*	trophon EPR, Sonex-HL*	trophon2, NanoNebulant*	trophon2, Sonex-HL*
5C1	✓	✓	✓	✓	✓	NA	NA	NA	NA
6C1	✓	✓	✓	✓	✓	✓	✓	✓	✓
7C2	✓	✓	✓	✓	✓	NA	NA	NA	NA
12L3	✓	NA	✓	✓	NA	NA	NA	NA	NA
11L4	✓	NA	✓	✓	✓	✓	✓	✓	✓
13L4	✓	✓	✓	✓	✓	✓	✓	✓	✓
14L4	✓	✓	✓	✓	✓	✓	✓	✓	✓
16L4	✓	✓	✓	✓	✓	✓	✓	✓	✓
18H5	✓	NA	✓	✓	NA	✓	✓	✓	✓
11M3	✓	✓	✓	NA	✓	✓	✓	✓	✓
L10-5v	✓	NA	✓	✓	✓	NA	NA	NA	NA
5P1	✓	✓	✓	✓	✓	✓	✓	✓	✓
8V4	✓	✓	✓	✓	✓	✓	✓	✓	✓
10V4	✓	✓	✓	✓	✓	✓	✓	✓	✓
P4-2	✓	✓	✓	✓	✓	✓	✓	✓	✓
9MC3	✓	✓	✓	✓	✓	✓	✓	✓	✓
10MC3	✓	✓	✓	✓	✓	✓	✓	✓	✓
10BCC3	✓	NA	✓	✓	NA	✓	✓	✓	✓
9VC2	✓	✓	✓	✓	✓	✓	✓	✓	✓
8VC3	✓	NA	NA	✓	NA	✓	✓	✓	✓
9VE3	✓	NA	NA	✓	✓	✓	✓	✓	✓
9VE4	✓	✓	✓	✓	✓	✓	✓	✓	✓
CW2	✓	✓	✓	✓	✓	NA	NA	✓	✓
CW5	✓	✓	✓	✓	✓	NA	NA	✓	✓
CW8	✓	✓	✓	✓	✓	NA	NA	NA	NA

✓ = approved

NA = not approved

An asterisk (*) indicates a reprocessing device with the approved disinfectant for use. For information about the reprocessing device, refer to the device manufacturer's instructions for use.

List of Approved Sterilization Systems and Sterilants

Follow the sterilization system manufacturer's instructions for sterilizing the transducer.

Note: An approved sterilant may discolor the transducer. There is no associated degradation of imaging performance or transducer reliability.

STERRAD 100NX Sterilizer System

- Compatible transducers: 18H5
- Manufacturer: Advanced Sterilization Products (Johnson & Johnson)
- Sterilization method/type: Hydrogen peroxide gas plasma
- Approved sterilization cycles:
 - Standard
 - Express

See also: For information about the sterilization cycles and required accessories, refer to the sterilization system manufacturer's instructions.

Transducer Accessories

⚠ WARNING: Ensure the accessories for a transducer are properly cleaned, disinfected, or sterilized as appropriate before each use to avoid possible patient contamination.

Attachment procedures for the following accessories are presented in this chapter or shipped separately with the device.

Accessory	Curved Array	Linear Array	Phased Array	<i>fourSight 4D</i>
Transducer Sheath	All	All	All	All
Biopsy Needle Guide Set	5C1	11L4	---	---
	6C1	12L3		
	7C2	13L4		
		14L4		
		L10-5v		
Needle Guides		16L4		
	9MC3	---	---	9VE3
	10MC3			9VE4
	10BCC3			

Transducer Sheaths (Transducer Covers)

-  **WARNING:** There have been reports of severe allergic reactions to medical devices containing latex (natural rubber). Health care professionals are advised to identify latex-sensitive patients and to be prepared to treat allergic reactions promptly. For additional information in the U.S.A., refer to FDA Medical Alert MDA91-1.
-  **WARNING:** To minimize the risk of cross-contamination and infectious diseases, a sterile, non-pyrogenic transducer sheath must be in place during procedures requiring sterility.
-  **WARNING:** Only a sterile transducer sheath provides the sterile barrier required for surgical or puncture procedures. To ensure sterility of a procedure, always place a sterile sheath on a transducer, as transducers cannot be sterilized using hot steam, cold gas, or Ethylene Oxide (EO) methods.
-  **WARNING:** After placing the sheath over the transducer, visually inspect the sheath to ensure there are no defects. Do not use the sheath if it has any holes or tears.
-  **Caution:** Siemens Healthineers recommends that you follow all instructions provided by manufacturers of sterile goods (transducer sheaths) to ensure proper handling, storage, and cycling of all sterile goods.

Transducer sheaths are single-use items used to ensure proper acoustic coupling and provide a prophylactic barrier for the intended ultrasound application. Sheaths are available for all transducers. Siemens Healthineers recommends the use of market-cleared transducer sheaths.


Using a disposable latex transducer sheath on a transducer reduces the possibility of cross-contamination. Always use a protective transducer sheath for endocavity exams, and when scanning an open wound or an area where the skin is not intact.

To place the transducer in a sheath:

Note: For sterile use, follow a sterile technique including use of sterile water and a sterile water-based coupling agent (gel), taking care not to contaminate the sheath.

Before applying any coupling agent (gel) to the sheath, remove any powder in the sheath by rinsing with water.

1. Remove the packaging and unfold the transducer sheath.
2. Apply a water-based coupling agent (gel) to the inside of the sheath and onto the face of the transducer.
3. Hold the transducer by the cable relief and unroll the sheath onto the transducer.
4. Pull the transducer sheath tightly over the face of the transducer to remove wrinkles.
5. Secure the sheath to the transducer housing or cable relief with the adhesive tapes or elastic bands.

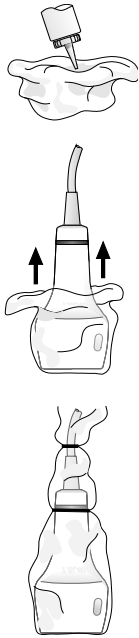

Example of placing a sheath over a transducer.

Disposal

While wearing protective gloves, remove the transducer sheath from the transducer. Dispose of the transducer sheath according to local, state, and regional laws and regulations for waste.

Storage

⚠ WARNING: Before use, examine sterile goods, such as sheaths, for any material flaws. Some packaging may list an expiration date. Any product showing flaws, or whose expiration date has passed, should not be used.

⚠ Caution: Do not store transducer sheaths in direct sunlight, as ultraviolet damage can result.

Latex products have a limited shelf life, and should be stored in a cool, dry, dark place with an ambient temperature between -5°C and +40°C and up to 80% relative humidity at +40°C.

Brackets and Needle Guides

⚠ WARNING: Percutaneous procedures always involve heightened risk to the patient and to the operator handling biopsy needle guides. Clinicians using Siemens Healthineers recommended biopsy devices under ultrasound guidance should be trained and must observe proper needle insertion sequencing with the needle guide in order to avoid undue discomfort and unnecessary risk and injury to the patient.

⚠ WARNING: If a reusable bracket becomes contaminated with tissue or fluids of a patient known to have Creutzfeldt-Jakob disease, then the bracket should be destroyed. Sterilization is not effective against Creutzfeldt-Jakob contamination.

⚠ WARNING: Reusable brackets are packaged non-sterile. Sterilize or high-level disinfect brackets products prior to their first use and after each use.

⚠ WARNING: The disposable needle guides are packaged sterile and are single-use items. Do not use if the packaging indicates signs of tampering or if the expiration date has passed.

Use the reusable bracket to attach needle guides to the transducer.

Refer to the in-box instructions for attachment and care procedures for all needle guide accessories, including information on reprocessing the reusable bracket.

Storage

Always clean and sterilize or high-level disinfect the reusable bracket after each needle puncture or biopsy procedure, before storing the reusable bracket.

7 Specialty Transducers


9VE3 Transducer	2
Imaging with the Transducer	2
9VE3 Technical Data	3
10BCC3 Biplane Transducer	4
Imaging with the Biplane Transducer	5
Activating Biopsy Guideline Function	5
10BCC3 Technical Data	6

9VE3 Transducer

⚠ WARNING: Before attempting to use endocavity transducers, you should be trained in ultrasonography and endocavity scanning techniques, and you should be thoroughly familiar with the safe operation of the ultrasound imaging system.

See also: For information about cleaning, disinfecting, and storing transducers, refer to Chapter 6 in this manual.

The 9VE3 curved array transducer is designed with *fourSight* 4D transducer technology.

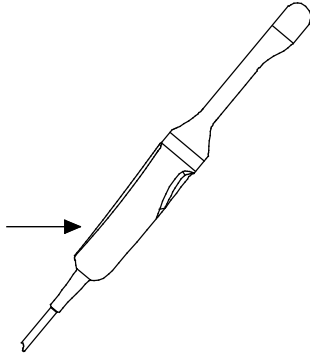

Example of the orientation indicator on the transducer handle.

Imaging with the Transducer

The transducer supports all the functions and controls that are available in 2D-mode, M-mode, Color, Power, Pulsed Doppler, and during volume acquisition.

To	Do This
Select a steering angle	Rotate [Steering Angle]. The steering angle increases or decreases relevant to the scan angle. An on-screen icon indicates the steering angle of the transducer when the steering angle is greater or less than 0°. <i>Example of icon during 2D-mode imaging.</i> <i>Example of icon during volume acquisition.</i>
Select a scan angle (to determine the amount of data for volume acquisition)	Rotate [Angle]. During 3-Scape imaging, the system computes the scan duration based on the selected angle, acquisition speed, and the current frame rate and displays the steering angle and acquisition time on the image screen.
Adjust the image quality	Rotate [Quality].

9VE3 Technical Data

See also: For information about operating frequencies, maximum display depth, environmental temperatures, and other technical descriptions, refer to Appendix A in this manual.


Transducer type:	Mechanically-driven curved array
Number of elements:	128
Focus:	35 mm $\pm 5\%$
Minimum depth:	3 cm
Maximum depth:	14 cm
Operating modes:	2D-mode, M-mode, 2D-mode with M-mode, 2D-mode with Doppler, Color, Power, Pulsed Wave Doppler
Maximum displayable field of view:	$145^{\circ} \pm 1^{\circ}$
Maximum scan angle (wobble angle):	$\pm 45^{\circ}$
Radius of curvature:	10.8 mm
Orientation:	On-screen graphical icon indicates steering angle

10BCC3 Biplane Transducer

⚠ WARNING: Before attempting to use endocavity transducers, you should be trained in ultrasonography and endocavity scanning techniques, and you should be thoroughly familiar with the safe operation of the ultrasound imaging system.

See also: For information about cleaning, disinfecting, and storing transducers, refer to Chapter 6 in this manual.

The 10BCC3 transducer is a biplane endocavity transducer for prostate exams.

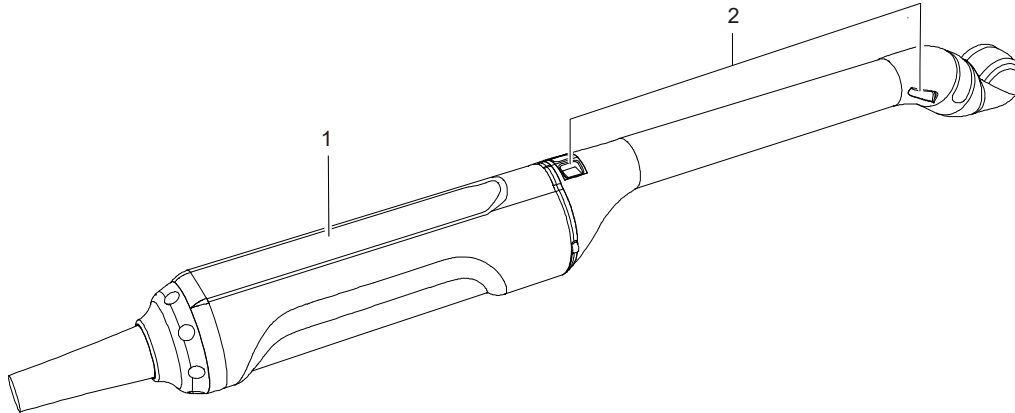

Example of indicators on the transducer.

- 1 Notches for indicating the transducer orientation
- 2 Ridge for attaching the needle guide

Imaging with the Biplane Transducer

Biplane imaging simultaneously scans two planes during real-time imaging. You can view either plane in a single display format or both planes on the image screen.

A scan plane icon indicates the transverse (T) and sagittal (S) planes.



Example of the transverse scan plane icon.



Example of the sagittal scan plane icon.

Use the configuration settings to specify the default display format for biplane and combined mode imaging.

To alternate between transverse or sagittal view in a single display format:

(Not available during Doppler imaging)

- Tap **Toggle View**. Or, press **UPDATE**.

To activate biplane in a dual display format:

- Tap **Live Bi-plane**.

The system activates a dual display format with sagittal and transverse views.

Note: Color mode is not available when **Live Bi-plane** is active.

To	Do This
Change the dual display format	▪ Tap Dual Layout and then select a setting. – Side by Side displays either real-time images or frozen images in a side-by-side format. – Top-Bottom displays either real-time images or frozen images in a vertical format.
Swap the transverse and sagittal views	▪ Press UPDATE. ▪ Tap Toggle View.
Show or hide a reference line at the intersection of the sagittal and transverse planes	▪ Tap Intersection Line.
Exit dual display format	▪ Tap Exit Dual.
Exit dual display format and activate real-time imaging	▪ Press 2D.

Activating Biopsy Guideline Function

See also: For information about the biopsy (puncture) guideline function, refer to Chapter 2 in the Clinical Applications Guide.

10BCC3 Technical Data

See also: For information about operating frequencies, maximum display depth, environmental temperatures, and other technical descriptions, refer to Appendix A in this manual.

Transducer type:	Curved array, biplane (sagittal and transverse)
Number of elements:	96 (each array)
Focus:	30 mm $\pm$ 5 mm
Minimum depth:	3 cm
Maximum depth:	14 cm
Operating modes:	2D-mode, Color, Color flow Doppler, M-mode, Pulsed Wave Doppler
Maximum field of view:	Sagittal plane: 164° Transverse plane: 164°
Radius of curvature:	10 mm $\pm$ 0.5 mm
Orientation:	On-screen graphical icons indicate sagittal and transverse scan planes

Appendix A Technical Description

Standard Features	3
Cybersecurity.....	3
System Operating Software.....	3
Monitor.....	3
Operator Control Panel.....	4
Mobility	4
Transducer Technology and Design Attributes.....	5
User-Accessible Connections.....	7
User-Accessible Cables	7
Gel Warmer	8
Storage Container	8
DICOM Connectivity	8
DICOM Structured Reporting.....	8
Left Ventricular Opacification (LVO)	8
Operating Modes.....	9
Formats	9
General Imaging Functions.....	10
2D-Mode Features	11
M-Mode Features	12
Pulsed Wave Doppler.....	12
Steerable Continuous Wave Doppler	13
Auxiliary Continuous Wave Doppler	13
Color Doppler	14
3-Scape Imaging	14
syngo eSieCalcs Software.....	15
fourSight 4D Ultrasound Imaging Technology	15
Speed of Sound.....	15
Clarify Vascular Enhancement (VE) Technology.....	15
Directional Power Doppler.....	15
Image Rotate during Pediatric Hip Measurement.....	15
Clips	16
Wireless.....	16
Storage and Archiving	16
Accessories and Options.....	17
List of Available Features	18
Options	20
Transducer Options.....	25
On-board Documentation Devices.....	27
Disposables	27
Consumables.....	27

Measurements, Calculations, and Reports	28
General Functions	28
2D-Mode Measurements and Calculations.....	28
M-Mode Measurements and Calculations	28
Doppler Measurements and Calculations.....	29
Exam-Specific Measurements and Calculations.....	29
Clinical Measurement Range and Accuracy	30
Clinical Measurements: Range and Accuracy	30
Virtual Touch Clinical Measurements: Range and Accuracy	34
Point Shear Wave Elastography (pSWE)	34
Ultrasound-Derived Fat Fraction (UDFF).....	34
System Requirements	35
Power Supply Requirements	35
Possible Combinations with Other Equipment.....	35
Leakage Currents	36
Input and Output Signal for Video and Data Transmission Connections	36
Video Standard.....	36
Wireless Network Connections.....	37
Environmental Requirements	38
Ultrasound System	38
Transducers	38
Protective Measures.....	39
Maximum Physical Dimensions.....	39
System Classifications.....	40
Standards Compliance.....	41
Radio and Telecommunications Standards	41

Standard Features

Cybersecurity

- Your ultrasound system includes a security package to protect patient confidentiality and system security.

See also: For information about cybersecurity, refer to Appendix C in this manual.

System Operating Software

- Windows-based operating software
- Multi-language capability, selectable during software installation
- User-modifiable configuration settings, including libraries of annotations and body markers
- Feature for defining specific exam, transducer, and image parameter settings
- Compatible with remote update handling for remote application support and remote troubleshooting

Monitor

- 54.5 cm (21.5 inch) high-resolution flat panel display (Monitor)
- 16:9 wide screen format
- 1920 x 1080 resolution
- Brightness controls
- $\pm 178^\circ$ wide viewing angle
- Energy saving feature enables a screen saver after a specified time period
- Variable positioning adjustments:
 - ACUSON Juniper system
 - Range of height: 141 cm to 163.4 cm; folds down for transport or portable exams to 122 cm
 - Left and right swivel: $\pm 90^\circ$
 - Tilt front ($+65^\circ$) and back (-15°)
 - ACUSON Juniper Select system
 - Range of height: 140 cm to 154 cm; folds down for transport or portable exams to 119 cm
 - Left and right swivel: $\pm 80^\circ$
 - Tilt front ($+80^\circ$) and back (-15°)
- Monitor arm adjustments:
 - ACUSON Juniper system
 - Left and right swivel: $\pm 320^\circ$
 - Vertical: up to 20°
 - ACUSON Juniper Select system
 - Left and right swivel: $\pm 80^\circ$
 - Horizontal: forward ($+30^\circ$) and backward (-30°)

Operator Control Panel

- Backlit for controls and keys
- Control panel height adjustment for standing and sitting positions
 - ACUSON Juniper system
 - Left/right swivel: $\pm 90^\circ$
 - Range of height: 11 cm
 - ACUSON Juniper Select system
 - Range of height: 11 cm
- Language-specific control panel
- Control panel layout supports ambidextrous operation
- Touch screen with brightness controls:
 - 33.8 cm (13.3 inch)
 - 1920 x 1080 resolution
 - Extended wide-angle viewing angle: $\pm 80^\circ$
 - Variable tilt angle: 40° to 60°
 - Language-specific virtual keyboard

Mobility

- ACUSON Juniper system
 - Four locking wheels
 - Directional swivel lock on each rear wheel
- ACUSON Juniper Select system
 - Four locking wheels
- Compact and lightweight industrial design
- Large steering handle for portability and maneuvering
- Transducer holders and cable management
- Tilt-down flat panel monitor
- Locking arm for flat-panel monitor

Transducer Technology and Design Attributes

- 1.5 to 14.5 MHz frequency range capability
- *fourSight* 4D transducer technology
- Wideband multiple frequency imaging
- Hanafy lens acoustic technology
- Virtual format imaging
- Transducer ports:
 - ACUSON Juniper system
 - Five ports for array transducers
 - One auxiliary continuous wave transducer port
 - ACUSON Juniper Select system
 - Four ports for array transducers
 - One auxiliary continuous wave transducer port

- Removable transducer holder for all transducer designs and gel bottle storage
- Adjustable temperature gel warmer connects to the ultrasound system's DC power socket and attaches to the side of the system in place of a transducer holder.
- Lightweight, ergonomic transducer design with SuppleFlex cables
- Range of operating frequencies

Transducer	2D-mode	Doppler
5C1	2.5 MHz - 5.0 MHz	2.0 MHz - 2.7 MHz
6C1	2.5 MHz - 5.0 MHz	2.0 MHz - 4.0 MHz
7C2	3.6 MHz - 5.7 MHz	2.7 MHz - 3.3 MHz
12L3	6.2 MHz - 13.3 MHz	4.0 MHz - 6.7 MHz
11L4	6.2 MHz - 10.0 MHz	4.0 MHz - 6.2 MHz
13L4	6.2 MHz - 13.3 MHz	4.0 MHz - 6.2 MHz
14L4	6.2 MHz - 14.0 MHz	4.0 MHz - 6.2 MHz
16L4	8.0 MHz - 13.3 MHz	5.3 MHz - 6.2 MHz
18H5	8.0 MHz - 14.5 MHz	6.2 MHz - 7.3 MHz
11M3	4.0 MHz - 8.0 MHz	3.6 MHz - 4.4 MHz
L10-5v	6.2 MHz - 10.0 MHz	4.0 MHz - 6.2 MHz
5P1	1.6 MHz - 4.0 MHz	1.8 MHz - 3.6 MHz
8V4	4.2 MHz - 7.2 MHz	3.6 MHz - 6.7 MHz
10V4	4.0 MHz - 8.9 MHz	4.0 MHz - 6.2 MHz
5VT	3.6 MHz - 5.7 MHz	2.9 MHz - 3.8 MHz
P4-2	2.5 MHz - 4.0 MHz	2.0 MHz - 3.3 MHz
9MC3	4.2 MHz - 7.3 MHz	4.0 MHz - 5.3 MHz
10MC3	5.0 MHz - 9.4 MHz	4.0 MHz - 6.2 MHz
10BCC3	5.7 MHz - 8.4 MHz	3.8 MHz - 5.0 MHz
9VC2	3.3 MHz - 6.1 MHz	2.7 MHz - 3.3 MHz
8VC3	3.3 MHz - 6.2 MHz	2.7 MHz - 3.3 MHz
9VE3	4.7 MHz - 7.3 MHz	4.0 MHz - 5.3 MHz
9VE4	4.7 MHz - 7.3 MHz	4.0 MHz - 5.3 MHz
CW2	---	2.0 MHz
CW5	---	5.0 MHz
CW8	---	8.0 MHz

- 2D-mode includes fundamental and harmonic imaging, without contrast
- Doppler includes pulsed wave, continuous wave, and color imaging

User-Accessible Connections

- On-board location for up to two image recording devices
- USB ports on the left side of the touch screen and back of the system for off-board printer connectivity, importing and exporting protocols, archiving data, and serviceability.
- Ethernet connection
- HDMI port for optional external monitor
- ECG connector socket
- Auxiliary ECG connector socket

User-Accessible Cables

Cable Type	Shielded	Length (m)	Ferrite
AC Mains Input (System power)	No	≤ 5.00 m	No
Footswitch	Yes	≤ 3.00 m	No
Barcode Reader	Yes	≤ 2.20 m	No
Gel Warmer	No	≤ 0.40 m	No
Black and White Printer (DC type)			
Power Cable (DC Input)	No	≤ 1.00 m	No
USB Cable	Yes	≤ 0.90 m	No
Color Printer			
Power Cable (AC Mains Input)	No	≤ 0.90 m	No
USB Cable	Yes	≤ 0.90 m	Yes
DVR			
USB Cable	Yes	≤ 0.90 m	Yes
Power Cable (AC Mains Input)	No	≤ 0.90 m	No
Power Cable (AC/DC Adapter)	No	≤ 1.85 m	Yes
HDMI Cable	Yes	≤ 1.20 m	Yes
ECG			
ECG Cable	Yes	≤ 2.95 m	Yes
ECG Leads	Yes	≤ 1.00 m	Yes
Transducer Cables			
5C1	Yes	≤ 2.10 m	Yes
6C1	Yes	≤ 2.20 m	Yes
7C2	Yes	≤ 2.10 m	Yes
12L3	Yes	≤ 2.10 m	Yes
11L4	Yes	≤ 2.10 m	Yes
13L4	Yes	≤ 2.20 m	Yes
14L4	Yes	≤ 2.20 m	Yes
16L4	Yes	≤ 2.10 m	Yes
18H5	Yes	≤ 2.10 m	Yes

Cable Type	Shielded	Length (m)	Ferrite
11M3	Yes	≤ 2.10 m	Yes
L10-5v	Yes	≤ 2.30 m	Yes
5P1	Yes	≤ 2.10 m	Yes
8V4	Yes	≤ 2.10 m	Yes
10V4	Yes	≤ 2.10 m	Yes
5VT	Yes	≤ 1.99 m	Yes
P4-2	Yes	≤ 2.10 m	Yes
9MC3	Yes	≤ 2.10 m	Yes
10MC3	Yes	≤ 2.10 m	Yes
10BCC3	Yes	≤ 2.30 m	Yes
9VC2	Yes	≤ 2.30 m	Yes
8VC3	Yes	≤ 2.10 m	Yes
9VE3	Yes	≤ 2.10 m	Yes
9VE4	Yes	≤ 2.30 m	Yes
CW2	Yes	≤ 1.83 m	Yes
CW5	Yes	≤ 2.10 m	Yes
CW8	Yes	≤ 2.20 m	Yes

Gel Warmer

- Heats the bottle of gel inserted into the gel warmer

Storage Container

- Basket attaches to the rear handle of the system.
- Side storage bin attaches to either the right or left side of the system

DICOM Connectivity

Note: The ultrasound system is the DICOM Storage Class User.

- Connectivity for digital data transfer to a DICOM network for storage
- Includes DICOM Worklist for querying and retrieving the patient worklist schedule from the hospital or radiology information system.

DICOM Structured Reporting

- Enables transfer of measurements to offline personal computers, workstations, and archiving systems
- Available only with the following exams: Cardiac, Obstetrics, Gynecology, Vascular

Left Ventricular Opacification (LVO)

- Provides improved resolution of endocardial borders and optimal evaluation for wall motion during contrast agent imaging
- Compatible transducers and exam types:
 - 5P1, P4-2: Cardiac

Operating Modes

See also: For a list of operating modes and combined modes, refer to Chapter 1 in this manual.

Formats

- Single and dual image display
- Display formats for M-mode and Doppler: side by side, 1/2 2D - 1/2 trace, 1/3 2D - 2/3 trace, 2/3 2D - 1/3 trace
- Full screen format for only the M-mode sweep or Doppler spectrum
- Virtual format imaging: linear, steered linear, sector, vector, or trapezoidal format
- Image display formats during review: 1x1, 1x2, 2x2, 3x3, 4x4, 5x5

General Imaging Functions

- User-adjustable transmit power: 0.20% to 100%
- User-adjustable field of view (FOV) format
- User-adjustable display depth: transducer dependent (1 cm to 35 cm)

Transducer	Maximum Display Depth (cm)
5C1	35
6C1	35
7C2	30
12L3	16
11L4	16
13L4	16
14L4	16
16L4	6
18H5	8
11M3	15
L10-5v	16
5P1	30
8V4	15
10V4	15
5VT	24
P4-2	30
9MC3	14
10MC3	14
10BCC3	14
9VC2	30
8VC3	30
9VE3	14
9VE4	14
CW2	---
CW5	---
CW8	---

2D-Mode Features

- Fundamental frequency, transducer dependent: up to 4
- Harmonic frequency, transducer dependent: up to 5
- Adjustable size and position of field of view
- Adjustable image size
- Gain in one decibel increments: -30 dB to 30 dB
- Dynamic range in one decibel increments: 10 dB to 90 dB
- Focal zones: up to 8
- Magnification: up to 10x
- Acquired 2D-mode frame rates, depending on the transducer and imaging depth: up to 2,070 fps (frames per second)
- CINE capture: up to 60 seconds
- Line density: 6 levels
- Persistence, transducer dependent: up to 6 levels
- Edge enhancement: 4 levels
- UltraArt universal image processing for speckle reduction, transducer dependent: up to 6 levels
- Grayscale maps, transducer dependent: up to 9
- Colorization maps: 16
- Compounding to enhance viewing of tissue differences by detection of subtle lesions and enhanced tissue differentiation
- Harmonics imaging to enhance visualization and reduce noise
- Auto Tissue Equalization (Auto TEQ) image optimization technology to optimize image brightness uniformity in the field of view by changing the DGC and overall gain
- Clarify vascular enhancement technology to increase the contrast resolution and enhance the boundary detection: 7 levels
- Speed of Sound to optimize the 2D image
- Wide Field of View (Wide FOV) to extend the field of view to 224°
- Left ventricular opacification (LVO) Contrast technology, which is a low power harmonic imaging intended for use with injectable ultrasound contrast agents to improve visualization of the endocardial borders and optimal evaluation for wall motion.

M-Mode Features

- Gain in one decibel increments: -30 dB to 30 dB
- Dynamic range in one decibel increments: 10 dB to 90 dB
- Multiple sweep speed selections: 10
- Edge enhancement settings: 4
- Fundamental frequency, transducer dependent: up to 4
- Harmonic frequency, transducer dependent: up to 5
- Grayscale maps, transducer dependent: up to 9
- Colorization maps: 16
- Anatomical M-mode for cardiac exams supports visualization of an M-mode sweep by angle or marker.

Pulsed Wave Doppler

- Gain in one decibel increments: -30 dB to 30 dB
- Dynamic range in one decibel increments, transducer dependent: 10 dB to 90 dB
- Simultaneous 2D-mode and Doppler display and 2D/Doppler with color (Triplex)
- Multiple sweep speed selections: 10
- Transducer frequencies, transducer dependent: up to 4
- Grayscale maps: 7
- Colorization maps: 12
- Angle correction in one degree increments: 0° to 89°
- Pulse Repetition Frequency (PRF) range: 152 Hz to 39,100 Hz
- Velocity scale range with 0 angle correction, transducer dependent: ± 0.8 cm/s to ± 840 cm/s
- Adjustable gate size: 0.5 mm to 20 mm
- Wall filter selections, depending on the transducer: up to 7
- Baseline shift: up to 13 levels
- Spectral invert
- System-generated trace function
- Spectral DTI (Doppler Tissue Imaging) to depict movement of tissue
- Auto Trace to trace the maximum velocity on a Doppler spectrum

Steerable Continuous Wave Doppler

(Available with the phased array transducer)

- Gain in one decibel increments: -30 dB to 30 dB
- Multiple sweep speed selections: 10
- Transducer frequencies: up to 3
- Grayscale maps: 7
- Colorization maps: 12
- Pulse Repetition Frequency (PRF) range: 152 Hz to 52,100 Hz
- Velocity scale range with 0 angle correction, transducer dependent: ± 1.15 cm/s to $\pm 1,100$ cm/s
- Wall filter selections, depending on the transducer: up to 7
- Baseline shift: up to 13 levels
- Spectral invert

Auxiliary Continuous Wave Doppler

- Gain in one decibel increments: -30 dB to 30 dB
- Multiple sweep speed selections: 10
- Transducer frequencies: 1
- Grayscale maps: 7
- Colorization maps: 12
- Pulse Repetition Frequency (PRF) range: 152 Hz to 52,100 Hz
- Velocity scale range with 0 angle correction, transducer dependent: ± 0.75 cm/s to $\pm 1,000$ cm/s
- Wall filter selections, depending on the transducer: up to 7
- Baseline shift: up to 13 levels
- Spectral invert

Color Doppler

- Gain in one decibel increments: -20 dB to 20 dB
- Pulse Repetition Frequency (PRF) range: 100 Hz to 25,500 Hz
- Velocity scale range, transducer dependent: ± 0.5 cm/s to 330.9 cm/s
- Adjustable color region of interest
- Levels of baseline shift: up to 13
- Transducer frequencies for color mode: up to 4
- Transducer frequencies for power mode: up to 4
- Color maps (color only): 10
- Power maps (power only), transducer dependent: up to 9
- Persistence: 5 levels
- Tissue/Priority: 5 selections
- Smoothing: 4 levels
- Wall filter: up to 4
- α -Blend: up to 4
- Color Doppler Tissue Imaging (DTI)
- Doppler Tissue Energy (DTE)

3-Scape Imaging

- The 3D Imaging feature acquires three-dimensional ultrasound images. Multiplanar rendering (MPR) provides a view of each segment of the volume as an arbitrary slice.
- Real-time reconstruction after free-hand acquisition
- Multiplanar rendering demonstrates imaging planes not accessible with normal scanning techniques

syngo eSie Calcs Software

- Provides area, circumference, and maximum diameter measurements using a border detection algorithm
- User defines the region of interest
- Can be used wherever manual trace is available

fourSight 4D Ultrasound Imaging Technology

- Compatible transducers and exam types:
 - 8VC3, 9VC2, 9VE3: Obstetrics, Early Obstetrics, Advanced Obstetrics, Abdominal, Gynecology, Pelvic Floor
 - 9VE4: Obstetrics, Early Obstetrics, Advanced Obstetrics, Gynecology, Pelvic Floor
- Provides general features for 4D imaging
- Provides advanced features for 4D imaging:
 - MultiSlice
 - Thick slice imaging (TSI)
 - Curved multiplanar rendering (MPR)

Speed of Sound

- Optimizes the 2D image by adjusting the speed of sound

Clarify Vascular Enhancement (VE) Technology

- Available in combination with 2D-mode, 2D-mode with Doppler, 2D-mode with Harmonic imaging, or compounding
- Uses flow information to decrease artifacts in the 2D-mode image, eliminating noise and reverberation from vessels, and resulting in an improved view of anatomical structures
- Decreased artifacts can enhance definition of both tissue and vessel walls by increasing contrast resolution and improving boundary detection

Directional Power Doppler

- Displays energy generated by the direction of blood flow

Image Rotate during Pediatric Hip Measurement

- Rotates the image acquired during pediatric hip angle measurement

Clips

- Prospective capture of motion image data
- Retrospective capture of motion image data
- Variable clip length (1 second to 60 seconds)
- AVI format
- DICOM format
- Supports storing of individual images within clips

Wireless

- Sends data, such as studies, images, and clips to a network location over a wireless network

Storage and Archiving

The ultrasound system supports storage and review of completed ultrasound studies, including static images, dynamic clips, measurements, calculations, and reports.

Accessories and Options

⚠ WARNING: Accessory equipment connected to the analog and digital interfaces must be certified according to the respective EN or IEC standards (for example, EN 60950 or IEC 60950 for data processing equipment and EN 60601-1 or IEC 60601-1 for medical equipment). Anyone who connects additional equipment to any of the signal input or signal output ports configures a medical system and is therefore responsible that the system complies with the requirements of the system standards EN 60601-1 or IEC 60601-1. Siemens Healthineers can only guarantee the performance and safety of the devices listed in the *Instructions for Use*. If in doubt, consult the Siemens Healthineers service department or your local Siemens Healthineers representative.

⚠ WARNING: You must only use the transducers, accessories, cables, and replacement parts for internal components specified by Siemens Healthineers to reduce the risk of increased RF (radio frequency) emissions or decreased immunity of the ultrasound system.

Note: Not all features and options described in this publication are available to all users. Please check with your Siemens Healthineers representative to determine the current availability of features and options.

Note: The system software with a power cord and one bottle of coupling agent (gel) are included with the ultrasound system.

The Siemens Healthineers-authorized accessories and options for your ultrasound system are listed in this chapter. The available options depend on the licenses purchased for your ultrasound system.

Note: To ensure compliance with the Medical Device Regulations, use only the devices listed in this chapter with your ultrasound system.

List of Available Features

The following table describes the availability of features for ACUSON Juniper and ACUSON Juniper Select ultrasound systems.

Feature	Ultrasound System	
	ACUSON Juniper	ACUSON Juniper Select
3-Scape Real-time Imaging	•	•
Auto IMT	•	✓
Battery Support	✓	✓
Clarify Vascular Enhancement (VE) Technology	•	•
Contrast Agent Imaging	✓	NA
Data Storage Encryption	✓	✓
DICOM Connectivity	•	•
DICOM Structured Reporting	•	•
Directional Power Doppler	•	•
DTI Doppler Tissue Imaging Capability	✓	✓
eSie Measure Workflow Acceleration Package	✓	✓
eSie Left Heart Measurement Package	✓	NA
<i>fourSight</i> 4D Ultrasound Imaging Technology	•	•
Image Rotate during Pediatric Hip Measurement	•	•
Left Ventricular Opacification (LVO)	•	•
LightSource	✓	✓
Needle Visualization	✓	NA
Panoramic Imaging	✓	✓
Protocols	✓	✓
QuikStart Standby Mode	✓	✓
Remote Services	✓	✓
Speed of Sound	•	•
Stress Echo	✓	NA
<i>syngo</i> Arterial Health Package	✓	✓
<i>syngo</i> Auto Follicle	✓	✓
<i>syngo</i> Auto OB Measurements	✓	✓
<i>syngo</i> eSie Calcs Software	•	•
<i>syngo</i> Velocity Vector Imaging Technology	✓	NA

Feature	Ultrasound System	
	ACUSON Juniper	ACUSON Juniper Select
Virtual Touch Point Shear Wave Elastography	✓	NA
Virtual Touch Strain Imaging	✓	NA
Virtual Touch UDFF (Ultrasound-Derived Fat Fraction)	✓	NA
Wireless	✓	✓
Alphanumeric Keyboard	✓	✓
Barcode Reader	✓	✓
Card Reader	✓	✓
ECG Leads	✓	✓
Footswitch	✓	✓
Gel Warmer	●	✓
Large Gel Holder	✓	✓
Physio Module	✓	✓
Storage Container	●	✓

● = included at purchase

✓ = available as an option at purchase

NA = not available

Options

Note: For information about the combination of options available for purchase, contact your Siemens Healthineers representative.

Options	Description
Panoramic Imaging	- Compatible with all imaging transducers (except the 10BCC3 transducer and CW and phased array transducers) - Panoramic imaging acquires two-dimensional ultrasound images with a composite, extended field of view
<i>syngo</i> Arterial Health Package*	(Requires the physio module option) - Compatible with 2D-mode clips or images acquired with a linear transducer - Provides a method of quantifying Carotid Intima-Media Thickness (CMT) - Uses semi-automated border detection to determine the maximum and average thickness of the intima-medial layer of the carotid artery
LightSource	(Requires the 3-Scape Imaging and <i>fourSight</i> 4D Imaging options) Provides detailed and near photorealistic views of the fetus
eSie Left Heart Measurement Package	(Requires the physio module option) - Compatible transducers and exam types: 5P1, P4-2: Cardiac - Clinical software program to perform left heart volume measurements using 2D-mode data - Traces the contours of the endocardial border of the left ventricle and the left atrium on the end-diastolic and end-systolic images in an apical four-chamber or apical two-chamber view - Includes methods for manually tracing the contours of the left ventricle or the left atrium - Displays the user-confirmed end-diastolic or end-systolic frames
Auto IMT	- Detects and measures the carotid intima-media thickness (IMT) in a user-positioned region of interest on a frozen 2D image acquired with a linear array transducer - Measures the maximum and average IMT of the carotid artery

Options	Description
<i>syngo</i> Auto Follicle	- Compatible transducers and exam type: 5C1, 6C1, 7C2, 9VC2, 8VC3, 10MC3, 9VE3, 9VE4, 9MC3: Gynecology - Provides system-detected distance measurements for user confirmation.
<i>syngo</i> Auto OB Measurements	- Compatible transducers and exam types: 5C1, 6C1, 7C2, 9VC2, 8VC3, 9MC3, 10MC3, 9VE3, 9VE4: Obstetrics, Early Obstetrics, Advanced Obstetrics 12L3, 11L4, 13L4, 14L4: Obstetrics - Provides system-detected measurements of specific fetal parameters for user confirmation.
<i>syngo</i> Velocity Vector Imaging Technology	- Compatible transducers and exam types: 5P1: Cardiac, Pediatric Echocardiography - P4-2: Cardiac 8V4, 10V4, 11M3: Neonatal Echocardiography, Pediatric Echocardiography 5VT: Transesophageal Echocardiography - Tracks and estimates tissue velocity and other motion and deformation parameters at selected points on a user-defined outline of a structure - Assists analysis of rotation, displacement, and radial strain of the left ventricle - Assists evaluation of pediatric or adult cardiac contraction by analyzing the systolic and diastolic ventricular strain and rotation
Stress Echo	(Requires the physio module option) - Compatible transducers and exam types: 5P1, P4-2: Cardiac
eSie Measure Workflow Acceleration Package	(Requires the physio module option) - Compatible transducers and exam types: 5P1: Cardiac, Pediatric Echocardiography - P4-2, CW2, CW5: Cardiac 8V4, 10V4, 11M3: Neonatal Echocardiography, Pediatric Echocardiography 5VT: Transesophageal Echocardiography - Reduces the keystrokes required to perform measurements in 2D-Mode, M-Mode, and pulsed wave Doppler
Protocols	- Provides a predefined checklist for a clinical workflow - Defines protocol steps (views) that include an image or clip with measurements, body markers, and annotations - Changes the sequence of views during a protocol - Pauses and resumes a protocol - Creates new workflow protocols or modifies existing protocols - Imports or exports protocols using a USB storage device

Options	Description
Virtual Touch Strain Imaging	- Provides a qualitative representation of relative tissue stiffness for the region of interest - Compatible transducers and exam types: 16L4, 12L3, 11L4, 13L4, 14L4, L10-5v: Breast, Thyroid 18H5: Thyroid
Virtual Touch Point Shear Wave Elastography	- Measures tissue shear velocity (Vs) and elasticity (E) for a selected region of interest using acoustic radiation force impulse (ARFI) technology of diagnostic ultrasound to induce tissue displacement. - Provides measurement labels for sites, masses, and liver segments - Compatible transducers and exam types: 5C1, 6C1: Abdominal, Abdominal Difficult
Virtual Touch UDFF (Ultrasound-Derived Fat Fraction)	- Measures tissue in a selected region of interest and calculates the Ultrasound-Derived Fat Fraction (UDFF) index for assessing steatotic liver disease - Measurements per label: up to 20 - Compatible transducers and exam types: 6C1: Abdominal, Abdominal Difficult
Needle Visualization	- Compatible transducers and exam types: 16L4: Bowel, Focused Assessment with Sonography for Trauma, Lung, Testis, Penile, Thyroid, Pediatric Hip Joint, Carotid, Breast, Musculoskeletal, Digital, Nerve 12L3, 11L4, 13L4, 14L4: Bowel, Focused Assessment with Sonography for Trauma, Lung, Obstetrics, Breast, Testis, Penile, Thyroid, Pediatric Hip Joint, Carotid, Peripheral Venous, Peripheral Arterials, Spine, Musculoskeletal, Digital, Nerve. - L10-5v: Bowel, Breast, Carotid, Digital, Lung, Musculoskeletal, Peripheral Arterials, Peripheral Venous, Testis, Thyroid 18H5: Testis, Penile, Thyroid, Musculoskeletal, Digital, Nerve, Intraoperative Abdomen, Carotid - Enhances the visualization of the needle during the biopsy or needle puncture procedures
Contrast Agent Imaging	Note: At the time of publication, the United States Food and Drug Administration has cleared ultrasound contrast agents only for use in left ventricular opacification to improve the delineation of the left ventricular endocardial border, liver for lesion characterization in adult and pediatric patients, and vesicoureteral reflux (VUR) in pediatric patients. Check the current regulation for the country where you are using this ultrasound system for contrast agent clearance. - Supports the following contrast agent destruction techniques: - Burst - Flash Sequencing - Compatible transducers: 6C1 14L4

Options	Description
DTI Doppler Tissue Imaging Capability	▪ Depicts tissue movement during pulsed wave Doppler mode
Battery Support	▪ Uses the ultrasound systems with battery power assistance. Allows a continuous scan without connecting the ultrasound system to the power supply outlet or even during power service interruptions.
ECG Leads	▪ ECG leads, standard U.S.A. ▪ ECG leads, standard European
Barcode Reader	▪ Provides direct entry of patient information in the registration form, for example, the patient name and ID, the performing physician, or the sonographer ▪ Connects to a USB port on the ultrasound system. ▪ Barcode Reader, DS8108
Card Reader	(Available only for systems in Japan) ▪ Provides direct entry of patient information in the registration form, for example, the patient name and ID, the performing physician, or the sonographer
Large Gel Holder	▪ Holder for storing large bottles of gel on the ultrasound system in place of a transducer holder
Alphanumeric Keyboard	▪ Retractable alphanumeric keyboard for attachment below the control panel
QuikStart Standby Mode	▪ For use during portable studies ▪ Decreases the time required to power on or off the ultrasound system using the installed battery to place the ultrasound system in a standby status
Footswitch	▪ Configurable dual watertight footswitch ▪ Requires connection to the ultrasound system with a USB connector ▪ Footswitch, FSU-2000U

Options	Description
Physio Module	▪ Physiologic function – ECG cable and leads – ECG auxiliary input/output cables – Stress adapter cable set – Auxiliary ECG function – Auxiliary signal – Respiratory trace ▪ Provides the ability to configure ECG capabilities for specialty applications that require display and monitoring of physiologic traces without continuous wave Doppler capabilities
Data Storage Encryption	 Caution: Data storage encryption on your ultrasound system requires a USB storage device with the recovery key for recovering encrypted data in the event of a system failure. The ultrasound system administrator must keep the recovery key in a secure location. ▪ Data storage encryption – Protects patient information and system settings stored on the ultrasound system by preventing unauthorized access – Includes a USB storage device with the recovery key for recovering encrypted data
Remote Services	▪ Enables a Siemens Healthineers representative to remotely view and control your system ▪ Requires using the phone number provided by your Siemens Healthineers representative to initiate remote assistance. The representative guides you through the setup of each communication session.

*Not available for sale from Siemens Healthineers.

Transducer Options

ACUSON Juniper System

Options	Description
Transducers, Curved Array	▪ 5C1 ▪ 6C1 ▪ 7C2 ▪ 11M3
Transducers, <i>fourSight</i> 4D	▪ 8VC3* ▪ 9VC2 ▪ 9VE3* ▪ 9VE4
Transducers, Linear Array	▪ 12L3 ▪ 11L4 ▪ 13L4* ▪ 14L4 ▪ L10-5v ▪ 16L4 ▪ 18H5
Transducers, Phased Array	▪ 5P1 ▪ 8V4 ▪ 10V4 ▪ 5VT ▪ P4-2
Transducers, Continuous Wave	▪ CW2 ▪ CW5 ▪ CW8*
Transducers, Endocavity	▪ 9MC3 ▪ 10MC3 ▪ 10BCC3
Transducer Accessories	▪ Transducer Sheaths (except CW transducers) ▪ Biopsy Needle Guide Set – Verza Guidance System: 6C1, 13L4, 14L4 – In-plane ultrasound needle guide, Ultra-Pro II: 5C1, 7C2, 11L4, 12L3, L10-5v – In-plane ultrasound needle guide, Infiniti Plus: 16L4 ▪ Reusable Endocavity Needle Guide: 9VE3, 9MC3, 10MC3, 10BCC3 ▪ Disposable Endocavity Needle Guide: 9MC3, 10MC3, 9VE4, 10BCC3

*Not available for sale from Siemens Healthineers

ACUSON Juniper Select System

Options	Description
Transducers, Curved Array	▪ 5C1
Transducers, <i>fourSight</i> 4D	▪ 8VC3*
	▪ 9VC2
Transducers, Linear Array	▪ 12L3
	▪ 11L4
	▪ L10-5v
Transducers, Phased Array	▪ 5P1
	▪ P4-2
Transducers, Continuous Wave	▪ CW2
	▪ CW5
Transducers, Endocavity	▪ 10MC3
Transducer Accessories	▪ Transducer Sheaths (except CW transducers) ▪ Biopsy Needle Guide Set – In-plane ultrasound needle guide, Ultra-Pro II: 5C1, 11L4, 12L3, L10-5v ▪ Reusable Endocavity Needle Guide: 10MC3 ▪ Disposable Endocavity Needle Guide: 10MC3

*Not available for sale from Siemens Healthineers

On-board Documentation Devices

Note: Up to two documentation devices (black and white printer and color printer or DVR) can be integrated into the system.

- Black and White Printer, P95DW-DC(S), Mitsubishi*
- Black and White Printer, UP-D898DC, Sony
- Black and White Printer, UP-D711, Sony
- Color Printer, CP30DW, Mitsubishi*
- Color Printer, UP-D-25MD, Sony
- DVR UR-50BD-SK Recorder (NTSC/PAL, 115V/230V), TEAC*

*Not available for sale from Siemens Healthineers

Disposables

- Disposable Transducer Sterile Cover
- Disposable ECG Electrodes
- 3-lead ECG Cable
- Stress Echo Adapter Cable Set
- Protective Overlay for Control Panel
- Disposable Needle Guides and Brackets

Consumables

- Contact Scan Gel
- Sterile Contact Gel
- Black and White Printer Paper

Measurements, Calculations, and Reports

Two categories of measurements and calculations are available on the ultrasound system: general and exam-specific.

General Functions

- 2D-mode, M-mode, and Doppler support up to twenty (20) measurements per image
- Label then Measure or Measure then Label workflows
- Customizable reports

2D-Mode Measurements and Calculations

- Distance
- Angle
- Circumference
- Area
- Volume
- Disk Volume
- Stenosis Ratio

M-Mode Measurements and Calculations

- Distance
- Time
- Heart Rate
- Slope
- Distance Ratio
- Time Ratio

Doppler Measurements and Calculations

- Velocity
- Heart Rate
- Time
- Time Ratio
- Acceleration
- Flow Volume
- Resistive Index (from peak systolic velocity and end diastolic velocity)
- Pulsatility Index (from system-generated trace or user-generated trace)
- Velocity Ratio
- Velocity Time Integral
- Pressure Half Time

Exam-Specific Measurements and Calculations

The measurement function is arranged by exam type and is available for use with all exam types. All exam types support the following measurement and report features:

- All general measurements and calculations
- Exam-specific patient report (editable)
- User-defined anatomy descriptions

See also: For information about exam-specific measurements and calculations, refer to the measurement chapters in the Clinical Applications Guide.

Clinical Measurement Range and Accuracy

The following table describes the variability in accuracy for each parameter:

Clinical Measurements: Range and Accuracy

The ultrasound system utilizes the selected speed of sound for all measurements. When the speed of sound is adjusted or an image preset is used with a different speed of sound, the anatomy measured at a different speed of sound produces different measured results. However, the following ranges and accuracy apply to all speeds of sound.

Direct Measurement	Range	Accuracy
Distance	0 cm – 35 cm	3% of the distance or 1.5 mm; whichever is greater assuming 1,540 m/sec speed of sound. Does not apply to trace tool. Distance tolerance using trace tool is user-dependent.
Distance using extended field of view	0 cm – 27.5 cm	Linear transducer: 5% of the distance or 2.5 mm; whichever is greater assuming 1,540 m/sec speed of sound. Curved transducer: 8% of the distance or 2.5 mm; whichever is greater assuming 1,540 m/sec speed of sound. Note: Measurement results on Panoramic images may vary depending on the user's scanning technique. Improper skill or technique can cause inaccurate measurements.
Distance using Auto IMT	Within the region of interest (ROI) box	±0.2 mm up to 10° perpendicular from the intima-media thickness (IMT) surface.
Trace Area	0 cm ² – 972 cm ²	6% of the area or 1.5 cm ² ; whichever is greater, assuming minimal operator error in tracing the desired object and assuming 1,540 m/sec speed of sound.
Time	0 sec – 18.0 sec	Less than ±1% of the sweep speed or ±10 msec, whichever is greater.
Velocity	0 cm/sec – 2200 cm/sec	10% of the velocity or 5 cm/s; whichever is greater using a calibrated flow phantom.
Physio		
Heart Rate	30 bpm – 300 bpm	Less than ±2 beats per minute or ±1%, whichever is greater.
Angle	0° – 180°	±2° of the intersecting line with 1 cm in minimum length.

2D-Mode Measurements

2D-mode Measurements	Range or Formula	Minimum Value	Maximum Value
Distance	0 cm – 35 cm	D-Tol(D)	D + Tol(D)
Distance Ratio	D1/D2	(D1-Tol(D1))/(D2+Tol(D2))	(D1+Tol(D1))/(D2-Tol(D2))
Mean Distance	(D1+D2+D3)/3	(D1-Tol(D1)+D2-Tol(D2)+D3-Tol(D3))/3	(D1+Tol(D1)+D2+Tol(D2)+D3+Tol(D3))/3
Ellipse Circumference	$\pi\{1/2[D1^2+D2^2]\}^{1/2}$	$\pi\{1/2[(D1-Tol(D1))^2+(D2-Tol(D2))^2]\}^{1/2}$	$\pi\{1/2[(D1+Tol(D1))^2+(D2+Tol(D2))^2]\}^{1/2}$
Ellipse Area	$\pi(D1)(D2)/4$	$\pi(D1-Tol(D1))(D2-Tol(D2))/4$	$\pi(D1+Tol(D1))(D2+Tol(D2))/4$
Trace Area	0 cm ² – 972 cm ²	A-Tol(A)	A+Tol(A)
% Stenosis by Area	100(A1-A2)/A1	100(A1-Tol(A1)-A2-Tol(A2))/(A1+Tol(A1))	100(A1+Tol(A1)-A2+Tol(A2))/(A1-Tol(A1))
% Stenosis by Diam.	100(D1-D2)/D1	100(D1-Tol(D1)-D2-Tol(D2))/(D1+Tol(D1))	100(D1+Tol(D1)-D2+Tol(D2))/(D1-Tol(D1))
Angle	0° – 180°	0°+Tol(Dg)	180°+Tol(Dg)
Strain Ratio	M _{ROI1} /M _{ROI2}	---	---

Legend

Abbreviation	Definition
D	Distance
Tol	Tolerance
A	Area
MROI	Mean strain in a region of interest
Dg	Degree

Doppler Measurements

Doppler Measurements	Range or Formula	Minimum Value	Maximum Value
Velocity	0 cm/sec – 2200 cm/sec	V-Tol(V)	V+Tol(V)
Time @ Slow Sweep	0 sec – 18.0 sec	T-Tol(T)@ Slow Sweep	T+Tol(T)@ Slow Sweep
Time @ Med Sweep	0 sec – 4.0 sec	T-Tol(T)@ Med Sweep	T+Tol(T)@ Med Sweep
Time @ Fast Sweep	0 sec – 1.0 sec	T-Tol(T)@ Fast Sweep	T+Tol(T)@ Fast Sweep
Time Averaged Maximum Velocity; TAMx	$(MxVel_1 + MxVel_2 + \dots + MxVel_N)/N$	$(MxVel_1 - Tol(MxVel_1) + MxVel_2 - Tol(MxVel_2) + \dots + MxVel_N - Tol(MxVel_N))/N$	$(MxVel_1 + Tol(MxVel_1) + MxVel_2 + Tol(MxVel_2) + \dots + MxVel_N + Tol(MxVel_N))/N$
Time Averaged Mean Velocity; TAMn	$(MnVel_1 + MnVel_2 + \dots + MnVel_N)/N$	$(MnVel_1 - Tol(MnVel_1) + MnVel_2 - Tol(MnVel_2) + \dots + MnVel_N - Tol(MnVel_N))/N$	$(MnVel_1 + Tol(MnVel_1) + MnVel_2 + Tol(MnVel_2) + \dots + MnVel_N + Tol(MnVel_N))/N$
Resistivity Index	(PS-ED)/PS	$(PS - Tol(PS) - ED - Tol(ED))/(PS + Tol(PS))$	$(PS + Tol(PS) - ED + Tol(ED))/(PS - Tol(PS))$
Pulsatility Index	(PS-ED)/TAV	$(PS - Tol(PS) - ED - Tol(ED))/(TAV + Tol(TAV))$	$(PS + Tol(PS) - ED - Tol(ED))/(TAV - Tol(TAV))$
Slope (Acceleration)	(V1-V2)/T	$(V1 - Tol(V1) - V2 - Tol(V2))/(T + Tol(T))$	$(V1 + Tol(V1) - V2 + Tol(V2))/(T - Tol(T))$
Systolic/Diastolic Ratio	S/D	$(PS - Tol(PS))/(ED + Tol(ED))$	$(PS + Tol(PS))/(ED - Tol(ED))$

Legend

Abbreviation	Definition
V	Velocity
sec	Seconds
Tol	Tolerance
T	Time
Slow Sweep	Slow sweep
Med Sweep	Medium sweep
Fast Sweep	Fast sweep
MxVel	Maximum velocity
N	Number of velocities averaged
MnVel	Mean velocity
MdVel	Mode velocity
PS	Peak systolic velocity
ED	End diastolic velocity
TAV	Time averaged velocity
S/D	Peak systolic velocity to end diastolic velocity calculation

M-Mode Measurements

M-Mode Measurements	Range or Formula	Minimum Value	Maximum Value
Distance	0 cm – 35 cm	D-Tol(D)	D+Tol(D)
Time @ Slow Sweep	0 sec – 18.0 sec	T-Tol(T) @ Slow Sweep	T+Tol(T) @ Slow Sweep
Time @ Med Sweep	0 sec – 4.0 sec	T-Tol(T) @ Med Sweep	T+Tol(T) @ Med Sweep
Time @ Fast Sweep	0 sec – 1.0 sec	T-Tol(T) @ Fast Sweep	T+Tol(T) @ Fast Sweep
Heart Rate	30 bpm – 300 bpm	N/Time for N beats - N/(Time for N beats+Tol(T)) N = 1, 2, 3, 4, or 5	N/Time for N beats - N/(Time for N beats-Tol(T)) N = 1, 2, 3, 4, or 5

Legend

Abbreviation	Definition
D	Distance
Tol	Tolerance
T	Time
Slow Sweep	Slow sweep
Med Sweep	Medium sweep
Fast Sweep	Fast sweep
N/Time for N beats	Number of heart beats/number of heart beats per minute
N = 1, 2, 3, 4, or 5	Number of heart beats = 1, 2, 3, 4, or 5

Virtual Touch Clinical Measurements: Range and Accuracy

Point Shear Wave Elastography (pSWE)

Transducer	Direct Measurement	Range	Tolerance Symbol	Tolerance Function or Value (Accuracy)
5C1	Shear Wave Velocity	1.0 m/s to 5.0 m/s	Tol(Vs)	±10% over the depth range of 2.0 - 8.0 cm
	Elasticity (Young's Modulus)	3.0 kPa to 75.0 kPa	Tol(E)	±20% over the depth range of 2.0 - 8.0 cm
6C1	Shear Wave Velocity	1.0 m/s to 5.0 m/s	Tol(Vs)	±10% over the depth range of 2.0 – 8.0 cm
	Elasticity (Young's Modulus)	3.0 kPa to 75.0 kPa	Tol(E)	±20% over the depth range of 2.0 – 8.0 cm

Ultrasound-Derived Fat Fraction (UDFF)

UDFF had strong performance as compared to magnetic resonance imaging proton density fat fraction (MRI-PDFF) with an area under the receiver operating characteristic (AUROC) of 0.85 or higher.

UDFF had good agreement with MRI-PDFF in adults*. The Bland-Altman analysis of UDFF and MRI-PDFF had limits of agreement between -4.2 to 10.0 percentage points, with a mean absolute difference of 2.95 percentage points.

Transducer	Direct Measurement	Tolerance Function or Value (Accuracy)
6C1	Ultrasound-Derived Fat Fraction (UDFF) index	The overall range values of attenuation coefficient (AC) and backscatter coefficient (BSC) values with UDFF are within ±20% of the nominal values of the phantoms over the depth range of 2.8 cm to 8.2 cm

*UDFF exhibits bias toward slightly larger values versus MRI-PDFF. Larger differences between UDFF and MRI-PDFF occur at higher UDFF and MRI-PDFF values, with some outliers present at these levels.

System Requirements

Power Supply Requirements

Requirement		Ultrasound System
Input	MAINS Voltage	100V~ to 240V~
	Maximum Current	5.9 amps
	Frequency	50/60 Hz
Accessory Outlet	Outlet Voltage	100V~ to 240V~
	Maximum Current	150 VA
	Frequency	50/60 Hz

Possible Combinations with Other Equipment

⚠ WARNING: Accessory equipment connected to the analog and digital interfaces must be certified according to the respective EN or IEC standards (for example, EN 60950 or IEC 60950 for data processing equipment and EN 60601-1 or IEC 60601-1 for medical equipment). Anyone who connects additional equipment to any of the signal input or signal output ports configures a medical system and is therefore responsible that the system complies with the requirements of the system standards EN 60601-1 or IEC 60601-1. Siemens Healthineers can only guarantee the performance and safety of the devices listed in the *Instructions for Use*. If in doubt, consult the Siemens Healthineers service department or your local Siemens Healthineers representative.

Only the peripheral devices listed in the *Instructions for Use* are approved for use with the ultrasound system. Any use of other devices with the system will be at the user's risk and may void the system warranty.

On-board peripheral devices must be installed by an authorized Siemens Healthineers representative or approved third party. Check with your sales representative.

Leakage Currents

⚠ WARNING: Connecting peripheral products and accessories from non-isolated sources may result in chassis leakage current exceeding safe levels.

Input and Output Signal for Video and Data Transmission Connections

Port	Location	Connector
HDMI	Input/Output panel	19-pin HDMI (1 output)
USB	Input/Output panel	Series A type (4 ports)
USB	Left side of touch screen	Series A type (2 ports)
Ethernet (10BaseT/100BaseT)	Input/Output panel	RJ45 (1 port)
Digital Video Interface (DVI)	Input/Output panel	24-pin DVI-type (1 output)
S-Video	Input/Output panel	4-pin MINI-DIN (1 output)

Video Standard

- HDMI
- Digital Video Interface (DVI)
- S-Video

Wireless Network Connections

The ultrasound system supports the following options for connection to wireless networks.

Network standard:	▪ 802.11a (5.150 GHz to 5.850 GHz) with Orthogonal Frequency Division Multiplexing (OFDM) at 12 dBm $\pm$2 dBm @ 54 Mbps ▪ 802.11b (2.400 GHz to 2.4835 GHz) with Direct Sequence Spread Spectrum (DSSS) at 15 dBm $\pm$2 dBm @ 11 Mbps ▪ 802.11ac (5.150 GHz to 5.850 GHz) with OFDM at 10 dBm $\pm$2 dBm @ MCS9 ▪ 802.11g (2.400 GHz to 2.4835 GHz) with OFDM at 13 dBm $\pm$2 dBm @ 54 Mbps ▪ 802.11an (5.150 GHz to 5.850 GHz) with OFDM at 11 dBm $\pm$2 dBm @ MCS7 ▪ 802.11gn (2.400 GHz to 2.4835 GHz) with OFDM at 12 dBm $\pm$2 dBm @ MCS7
Authentication:	▪ WPA ▪ WPA2 ▪ WPA PSK ▪ WPA2 PSK ▪ Open
Encryption:	▪ None ▪ TKIP ▪ AES
Extensible Authentication Protocol (EAP):	▪ PEAPv0 (PEAP-MSHCAPv2) ▪ TLS

Environmental Requirements

Electromagnetic Compatibility (EMC) Note: Operating the ultrasound imaging system in close proximity to sources of strong electromagnetic fields, such as radio transmitter stations or similar installations, as well as portable and mobile RF (radiated frequency) may lead to temporary degradation or interference visible on the monitor screen. A lightening of image background may be noticed while visualizing hypoechoic structures, or color spectral interference, or jitter, or horizontal lines in the image screen may occur. However, the device has been designed and tested to withstand such interference and will not be permanently damaged. Wait until the interference disappears or restart the ultrasound system.

Ultrasound System

	During Operation	During Storage or Transportation
Atmospheric pressure:	700 hPa to 1060 hPa	500 hPa to 1060 hPa
Relative humidity:	10% to 80%, non-condensing	10% to 95%, non-condensing
Temperature:	---	---
System without a battery pack	+10°C to +40°C	-20°C to +60°C
System with a battery pack	+10°C to +35°C	-20°C to +50°C

Transducers

	During Operation	During Storage or Transportation
Atmospheric pressure:	---	---
5C1, 6C1, 7C2, 12L3, 11L4, 13L4, 14L4, 16L4, 18H5, L10-5v, 11M3, 8V4, 10V4, 5VT, P4-2, 9MC3, 10MC3, 10BCC3, 8VC3, 9VE3, CW8	700 hPa to 1060 hPa	500 hPa to 1060 hPa
5P1, 9VC2, 9VE4	500 hPa to 1060 hPa	500 hPa to 1060 hPa
CW2, CW5	---	680 hPa to 1060 hPa
Relative humidity:	---	---
5C1, 6C1, 7C2, 12L3, 11L4, 13L4, 14L4, 16L4, 18H5, L10-5v, 11M3, 8V4, 10V4, 5VT, P4-2, 9MC3, 10MC3, 10BCC3, 8VC3, 9VE3, CW8	10% to 80%, non-condensing	10% to 95%, non-condensing
5P1	Less than 95%, non-condensing	Less than 95%, non-condensing
9VC2, 9VE4	20% to 80%, non-condensing	20% to 80%, non-condensing
CW2, CW5	---	Less than 80%, non-condensing
Temperature:	---	---
5C1, 6C1, 7C2, 12L3, 11L4, 13L4, 14L4, 16L4, 18H5, L10-5v, 8V4, 10V4, P4-2, 10MC3, 10BCC3, CW2, CW5, CW8	+10°C to +40°C	-40°C to +60°C
5P1, 5VT	+10°C to +40°C	-10°C to +50°C
8VC3, 9VE3	+10°C to +40°C	-30°C to +70°C
9VC2, 9VE4	+20°C to +40°C	-10°C to +50°C
9MC3, 11M3	+10°C to +40°C	-20°C to +60°C

Note: Needle guides are excluded from the environmental requirements. Refer to the ranges included on the manufacturer's label.

Protective Measures

Explosion protection: This product is not designed for operation in areas subject to explosion hazards.

Maximum Physical Dimensions

ACUSON Juniper System

Width: 585.3 mm
Height: 1634 mm
Depth: 845.1 mm
Weight: 76 kg without documentation devices
Not to exceed 100 kg

ACUSON Juniper Select System

Width: 585.3 mm
Height: 1540 mm
Depth: 845.1 mm
Weight: 70 kg without documentation devices
Not to exceed 100 kg

System Classifications

- Type of protection against electrical shock:
 - Class I, external powered
- Degree of protection against electrical shock:
 - Type BF applied part for transducers
 - Type CF defibrillation-proof applied part for ECG
- Degree of protection against harmful ingress of water:
Ordinary equipment (transducers and footswitch: IPX8)
- Degree of safety of application in the presence of a flammable anesthetic material with air or with oxygen or nitrous oxide:
Equipment not suitable for use in the presence of a flammable anesthetic mixture with air or with oxygen or nitrous oxide.
- Mode of operation:
Continuous operation

Standards Compliance

The diagnostic ultrasound system is in compliance with the following standards, including all applicable amendments at the time of product release.

Quality Standards

- FDA QSR 21 CFR Part 820
- EN ISO 13485 and ISO 13485
- ISO 9001

Design Standards

- ANSI/AAMI ES 60601-1
- CAN/CSA-C22.2 No. 60601-1
- EN 60601-1 and IEC 60601-1
- EN 60601-1-2 and IEC 60601-1-2 (Class B)

Note: The system is a Class A device when the barcode reader, DVR, or 9VC2, 10BCC3, or 9VE4 transducer is in use.

- EN 60601-1-6 and IEC 60601-1-6
- EN 60601-2-18 and IEC 60601-2-18
- EN 60601-2-37 and IEC 60601-2-37
- EN 62304 and IEC 62304
- EN 62366-1 and IEC 62366-1
- EN ISO 14971 and ISO 14971

Acoustic Output Standards

- IEC 62359, Test Methods for the Determination of Thermal and Mechanical Indices Related to Medical Diagnostic Ultrasonic Fields
- AIUM/NEMA UD-2, Acoustic Output Measurement Standard for Diagnostic Ultrasound Equipment
- AIUM/NEMA UD-3, Standard for Real-Time Display of Thermal and Mechanical Acoustic Output Indices on Diagnostic Ultrasound Equipment

Radio and Telecommunications Standards

- CFR 47 FCC Part 15.247
- CFR 47 FCC Part 15.107
- CFR 47 FCC Part 15.109
- ETSI EN 300 328
- ETSI EN 301 489-1
- ETSI EN 301 489-17
- ETSI EN 301 893

Appendix B Patient Data

Patient Registration Form	2
Patient Information	2
Exam	3
Patient History	3
Institution	4
Patient Browser	5
Study List.....	6
Searching for a Patient Study.....	6
Querying and Retrieving Patient Data	7
Selections for Exporting Patient Studies.....	8
Selections for Managing Transfer Jobs	9
Worklist	10
Searching for a Patient on a Worklist Server.....	10
Editing MPPS Data	11
Appending a Scheduled Procedure Step.....	11
Adjusting the Study List	12

Patient Registration Form

The patient registration form contains general and exam-specific information. Use the configuration settings to select a long version of the patient registration form to include all fields on the form or to select a short version of the form.

Note: The patient name and ID can include letters, numbers, special characters, and single quotes (').

Patient Information

Field	Description
Last Name	Patient's last name
First Name	Patient's first name
Middle Name	Patient's middle initial(s)
Pre	List of prefixes, for example, Dr., Mrs., Ms.
Suf	List of suffixes, for example, Jr., II
Patient ID	Identification code for the patient. If you do not enter a patient identification code, the system generates a unique identifier starting with the system date and time. This code displays on the image screen. If there is not enough space for the entire code to display on the image screen, the system displays part of the code followed by three dots (...).
Date of Birth	Patient's date of birth using the date format specified in the configuration settings. When you enter or change this value, the system calculates and displays the value for Age .
Age	The system calculates and displays the patient's age based on the date of birth, if entered.
(Gender)	Patient's gender. If no gender is selected, the ultrasound system assigns gender based on the exam. ▪ Female for the obstetric, early obstetric, and gynecological exams ▪ Male is not assigned for any exams ▪ Other for all exams other than obstetric, early obstetric, and gynecological exams
Height	Patient's height
Weight	Patient's weight
BSA	Patient Body Surface Area in m². The system calculates the body surface area (BSA) using metric or U.S. formulas based on height and weight, as specified in configuration settings. ▪ Metric: $BSA = 0.007184 \times (\text{weight})^{0.425} \times (\text{height})^{0.725}$ ▪ U.S. $BSA = 0.007184 \times (\text{weight} \times 0.454)^{0.425} \times (\text{height} \times 2.54)^{0.725}$ The BSA is included in the patient report.

Exam

Field	Description
Transducer	List of the available transducers
Exam	List of the exam types available for the selected transducer
Protocol	List of workflow protocols
Requested Procedure	Requested procedure from the worklist server
Scheduled Procedure(s)	Scheduled procedure step from the worklist server
Request ID	Identification code for the procedure in the imaging service request. An imaging service request identifies attributes that are common among the requested procedures. Typically generated by a Hospital Information System/Radiology Information System server.
Accession No.	Identification code used for billing purposes indicating the sequence of the current study as related to other studies for this patient. Typically generated by a worklist server, for example, a Hospital Information System/Radiology Information System server.
Indication	Information describing the symptom or circumstance regarding the advisability or necessity of a specific medical procedure

Patient History

Field	Description
General	General
BP	Patient's blood pressure in mmHg
Additional Information	Additional information about the patient
OB	Obstetric
Date	Start date of the patient's last menstrual period or the date of in-vitro fertilization
LMP	Last menstrual period Double asterisks (**) before the estimated value indicate a calculated value using the estimated date of delivery.
EDD	Estimated date of delivery Calculation of the date is based on entries for last menstrual period or date of in-vitro fertilization.
Gravida	Patient's history of pregnancies
Aborta	Enter the number for each selection.
Para	These entries are transferred to the patient report but do not display on the image screen.
Ectopic	
Fetal Age	Calculation of the estimated fetal age in weeks and days is based on entries for last menstrual period, date of conception, or estimated date of delivery.
Additional Information	Additional information about the patient

Field	Description
GYN	Gynecology
Date	Start date of the patient's last menstrual period or the date of in-vitro fertilization
Gravida	Patient's history of pregnancies
Aborta	Enter the number for each selection.
Para	These entries are transferred to the patient report, but do not display on the image screen.
Ectopic	
Additional Information	Additional information about the patient
Vascular	Vascular
BP (Right)	Blood pressure of the right arm using the systole over diastole measurement in mmHg
BP (Left)	Blood pressure of the left arm using the systole over diastole measurement in mmHg
ABI (Left/Right)	Ankle Brachial Index
Additional Information	Additional information about the patient
Urology	Urology
PSA	Prostate Specific Antigen
Additional Information	Additional information about the patient


Institution

Use the configuration settings to update the lists of physicians and operators.

Field	Description
Performing Physician	Identifying information for the performing physician
Referring Physician	Identifying information for the referring physician
Operator	Identifying information for the operator

Patient Browser

Use the patient browser to search for patient records and view patient data. The patient browser displays a list of patient studies in the selected storage location, for example, the local database (hard disk), a configured server, or a USB storage device. You can also delete studies from the patient browser and export or copy studies to external media.

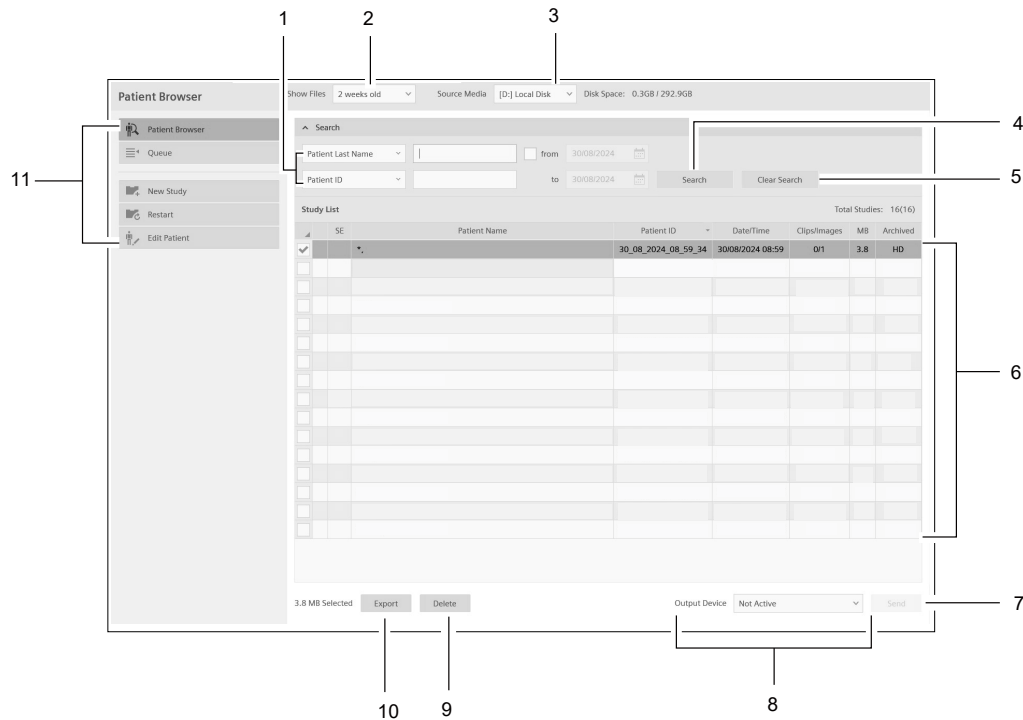

Example of the patient browser.

- 1 **Search criteria**
Use these fields to filter the list of studies. For example, you can search by patient name, study type, or a range of dates.
- 2 **Show Files**
Select the date range of studies to include in the search. For example, you can display all studies or only studies performed within the last month.
- 3 **Source Media**
Select the storage location of studies to include in the search.
- 4 **Search**
Click this button to activate the search based on the search criteria.
- 5 **Clear Search**
Click this button to display all studies stored in the selected location and to remove the search criteria.

6 Study List

Use this list to view study information, select a study for review, or restart a study. The system displays a graphical indicator to the left of the study currently open.

7 Send

Click this button to export selected studies to the selected external device or printer. For example, you can send selected studies to a DICOM storage server, a network export server, or a printer.

8 Output Device

Select the export location for the selected studies. For example, you can send selected studies to a DICOM storage server, a network export server or a printer.

9 Delete

Click this button to delete the selected studies from the system hard disk.

10 Export, Import

Click this button to display selections for exporting or importing studies.

11 Navigation menu

Use these controls to display other patient information screens, for example, the DICOM queue or the patient registration form.

Study List

Column Name	Description
Patient Name	Patient name entered on the patient registration form
Patient ID	Patient identifier entered on the patient registration form
Date/Time	Date and time the study was performed
Clips/Images	Number of clips and images saved with the study. For example, 4/18 indicates the study includes 4 clips and 18 static images.
MB	Size of the study in megabytes
Archived	Location of the archive of saved studies, for example, HD (hard disk) or USB
SE	Status of the stress echo protocol, for example, active, paused, or completed.

Searching for a Patient Study

You can search for a patient study saved on the local database and connected external storage devices.

To find a patient study:

Prerequisite: Access the patient browser.

1. Click **Source Media** and then select the required storage device.
2. Click **Show Files** and then, if necessary, specify the date range of the previous study.
3. Enter the required search criteria.
The system filters the list of patients as you enter the search criteria.
4. To clear the selected search criteria, click **Clear Search**.

Querying and Retrieving Patient Data

You can query and retrieve patient data from a connected DICOM archive using DICOM Query/Retrieve service.

Use the configuration settings to connect a storage server to the ultrasound system.

To retrieve a patient study from the storage server:

Prerequisite: Access the patient browser.

1. Tap **Query/Retrieve**.
2. Select a server from **Server List**.
3. Select the search criteria.

To	Do This
Confirm the connection to the network server	▪ Tap Ping.
Confirm the connection to the DICOM storage server	▪ Click Echo.
Specify detailed search criteria	1. Select the search criteria, for example, patient name or identification number. 2. Enter the required information in the search fields.
Select a date or range of dates for the patient studies	1. Click the check box next to the date input field. 2. Specify the date range for the studies. – Enter a start date and end date in the date input fields. – Click the calendar icons and then select a start date and end date.
Search for patient studies	▪ Select the required check box: – US Only searches all ultrasound systems connected to the storage server. – This system only searches only this ultrasound system.
Query the selected database	▪ Click Query.

4. Select the required study from the study list.
5. Select the required series from the series list.
6. Tap **Retrieve**.

Selections for Exporting Patient Studies

Selection	Description
Destination	Identifies target location for the exported studies
Identification	Identification
De-identify	Removes patient identifying information from the patient study Note: You can enter text to replace the patient name in the text box.
Format	Format
DICOM	Exports images in the DICOM (Digital Imaging and Communications in Medicine) format
Tiff/Avi	Exports images in Tagged Image File Format and exports clips in Audio Video Interleaved format
Option	Option
Finalize	Prevents additional storage to the inserted media
Eject	Ejects the media from the external device after the export is complete when the check box is selected
Eject	Ejects the media from the external device when the button is pressed
Status	Displays export progress
Export	Confirms the export is complete
Close	Closes the export dialog box
Send	Sends the study to the DICOM server

Selections for Managing Transfer Jobs

Use the queue in the patient browser to view and manage studies including images sent to a DICOM storage and a DICOM print server. Page tabs display along the top of the DICOM Store Queue and DICOM Print Preview.

Selection	Description
DICOM Store Queue	DICOM Store Queue
Type	Operation type, for example, Store (storage) or Commit (storage commitment)
Patient Name	Patient name entered on the patient registration form
Server	Destination server to receive the study
Time Sent	Date and time when the job was sent to the storage server
Information	Size of the files and the number of images or studies
State	Archive status of the study, for example, None (not archived) or Stored (archived)
Status	Information about the job, for example, in-progress or successful
Details	Error description if the job failed Note: SR indicates the content is a structured report.
Refresh Job	Updates the status information for all jobs in the queue
Retry Job	Resends the selected jobs, for example, if the transfer failed
Delete	Removes the selected jobs from the queue when the button is pressed
Delete All	Removes all jobs from the queue when the button is pressed
DICOM Print Preview	DICOM print preview
Printer Alias	Text entry to identify the configured device
Printer Type	Print output, for example, monochrome or color
Film size	Size of the print output
Orientation	Print orientation, for example, portrait or landscape
Copies	Text entry for the number of copies to print
Layout	Layout for the print output
Delete	Removes the image outlined from the layout page.
Cut	Outlines the selected image in yellow.
Paste	Inserts the image.
Print Page	Sends the current page of images to the queue
Print All	Sends all pages of images to the queue

Worklist

During patient registration, you can access patient studies from a worklist server.

Use the configuration settings to connect a worklist server to the ultrasound system.

Searching for a Patient on a Worklist Server

(Available when a worklist server is connected to the ultrasound system)

To search for a patient study:

Prerequisite: Access the worklist screen.

1. Select a server from the **Server** list.
2. To refresh the list of patients, click the refresh button.



3. Enter search criteria in the **Search** box or patient information fields.
The system filters the list of patients as you enter the search criteria.

To perform an advanced search for a patient study:

Prerequisite: Access the worklist screen.

1. Click **Advanced Search**.
2. Select a server from the **Server** list.
3. Select the search criteria.

To	Do This
Confirm the connection to the network server	▪ Click Ping.
Confirm the connection to the worklist server	▪ Click Echo.
Select a date or range of dates for the scheduled procedures	▪ Enter dates in the Date Range fields.
Search for scheduled procedures	▪ Select the required check box: – US Only searches all ultrasound systems connected to the worklist server. – This system only searches only this ultrasound system.
Enter detailed search criteria	▪ Enter the required information in the search fields. The system filters the list of patients as you enter the search criteria.
Clear all entries and selections	▪ Click Clear All.

4. Click **Get Worklist**.

Editing MPPS Data

You can edit the Modality Performed Procedure Step (MPPS) data, for example, descriptions and action codes.

Use the configuration settings to define and enable an MPPS server.

To edit MPPS data:

Prerequisite: Connect the ultrasound system to an MPPS server.

1. Press **PATIENT**.
2. Tap **MPPS**.
3. Enter or modify the required field, code meaning, or meaning.
4. Tap **Save**.

Appending a Scheduled Procedure Step

You can append a scheduled procedure step to a previous study.

Use the configuration settings to define and enable an MPPS server.

To append a procedure to an existing study:

Prerequisite: Configure an MPPS server.

1. Press **PATIENT**.
2. Tap **MPPS**.
3. To send the progress of the Performed Procedure Step (PPS) data:
 - a. Click **Click here to append**.
 - b. Enter the required field, code meaning or meaning, if necessary.
 - c. Click **Append**.
 - d. Click the required option.
 - **Completed** sends a successfully completed status to the DICOM MPPS server for the current procedure.
 - **Discontinued** sends a cancelled status to the DICOM MPPS server for the current procedure.
4. To exit without making any changes, tap **Cancel**.
The ultrasound system activates real-time imaging.

Adjusting the Study List

The study list displays the list of patient data for the current location, for example, the local database or a connected server.

You can adjust the study list by sorting the contents of the list, adding and removing columns, and changing the column widths.

To adjust the study list:

Prerequisite: Access patient registration. Access the worklist screen.

To	Do This
Sort the study list	Click a column heading.
Hide the list of patient studies	Click the hide button. 
Show the list of patient studies	Click the show button. 
Change the width of a column	- Access the pointer and then roll the trackball to position the pointer over the required column heading until the resize icon displays in the column.  - Press a trackball key and then drag the icon to resize the column.
Change the order of column headers	Click the required column heading and drag it to a new position.
Restore the study list to the default settings	- Click the additional selections button.  - Click Restore Defaults.

Appendix C Cybersecurity

Overview.....	2
Ultrasound System Security	2
Anti-Virus Protection.....	2
Firewall Policy	2
Protecting Patient Confidentiality.....	3
Security Responsibilities and Guidance	4
Recommendations for Your Hospital, Clinic, or Facility	4
Preventing, Detecting, and Removing Malware.....	5
Hardening the Ultrasound System.....	6
Managing User Accounts and Configuring Roles	7
Accessing the Ultrasound System for Emergency Use	8
Preventing Physical Access	8
Removing or Hiding Protected Health Information	9
Assuring Integrity and Confidentiality during Data Transmission.....	9
Backing Up Data and Disaster Recovery	10
Configuring Audit Controls and Managing Log Files	11
Third-Party Components.....	12
Remote Services.....	12

Overview

At the time of publication, the following information and recommendations provide measures for helping secure the operating environment of the ultrasound system and prevent unauthorized access.

Ultrasound System Security

The ultrasound system protects against security threats without disrupting the workflow or performance of the ultrasound system.

Anti-Virus Protection

The ultrasound system uses allowlisting software as an anti-virus solution. The anti-virus protection is preinstalled and activated on the ultrasound system before shipment. Allowlisting provides the following protection against security threats and vulnerabilities.

- Prevents execution of all unauthorized software and scripts
- Prevents modification to dynamic-link libraries (DLLs)
- Runs only trusted applications on the ultrasound system
- Accepts ultrasound system product updates through an authorized process
- Protects against threats from networks and infected USB storage devices

See also: Preventing, Detecting, and Removing Malware, page C-5

Firewall Policy

The ultrasound system is configured to use the operating system firewall as an added layer of protection.

Use the configuration settings to open or close communication through a port on the ultrasound system.

See also: For information about configuring the ports, refer to the security checklist in Chapter 1 in the System Reference.

Protecting Patient Confidentiality

The ultrasound system integrates security capabilities and security controls for assuring confidentiality of protected health information (PHI) before, during, and after an exam.

Patient data on the ultrasound system includes specific private, patient-identifying information, for example, patient demographics, biometric data, and information used for diagnosis, including text comments entered on images captured on the system and in patient reports.

- The image screen displays the patient data.
- Users can generate a hardcopy and electronic copy of images and reports.
- Users can also retrieve patient data from or save to removable media or a device on the network.

The ultrasound system can temporarily maintain the privacy of patient data in volatile memory until the system is powered off. Patient data can be persistently stored on local media and imported or exported through wired or wireless networks or through connected storage devices. The ultrasound system maintains the privacy of patient data during power service interruptions.

- You can protect patient confidentiality by using the security capabilities described in this appendix.

Security Responsibilities and Guidance

According to international standards for medical devices, the security of a medical device is the responsibility of the medical device manufacturer and your hospital, clinic, or facility as the delivery organization for healthcare.

- The ultrasound system has an underlying operating system. The operating system may be vulnerable to typical attacks, such as malware or denial-of-service.
- The following responsibilities of the medical device manufacturer prevent installation of unapproved software on the ultrasound system.
 - At the time of shipment to your hospital, clinic, or facility, the ultrasound system is current with applicable security patches.
 - Security updates are periodically offered for the ultrasound system to address critical vulnerabilities.

Use the configuration settings on the ultrasound system to view the Windows hotfixes installed on the system.

See also: For additional information, refer to the security checklist in Chapter 1 in the System Reference.

 - At the time of initial installation of the ultrasound system in your hospital, clinic, or facility, your service representative will confirm the integrity of the ultrasound system is the same as documented prior to shipment from the medical device manufacturer.
- Unintended consequences of data loss, misuse, or corruption not under the control of the ultrasound system, for example, after electronic communication from the ultrasound system to an information network or to storage devices are the responsibility of your hospital, clinic, or facility.
- You are responsible for the physical protection of the ultrasound system, for example, operating the system in a room with access control. Physical access to the system by unauthorized users permits the execution of exploits. The system contains patient data and should be protected against tampering and theft.
- Do not use removable media containing executable files from unknown sources.
- Do not load unnecessary external files onto the ultrasound system.

Recommendations for Your Hospital, Clinic, or Facility

- Apply appropriate administrative and physical safeguards to ensure the physical and electronic security of the network in your hospital, clinic, or facility, and with all medical devices on the network.
- Restrict access to medical devices by authorized users.
- Apply a firewall between the network of your hospital, clinic, or facility and external networks. Optionally, connect a hardware firewall between the network of your hospital, clinic, or facility and the ultrasound system.
- Monitor security incidents on your network.
- Prevent access of medical devices to the internet.
- Read the topics on the following pages of this appendix for an overview of the security capabilities of the ultrasound system. Each topic contains additional recommended actions for your hospital, clinic, or facility.

Preventing, Detecting, and Removing Malware

The ultrasound system supports an anti-malware protection mechanism for the prevention, detection, and removal of malware.

- The ultrasound system incorporates McAfee Application Control; only software authorized by the medical device manufacturer is executable.
- The ultrasound system captures security-related activities in a security log file.
- Ultrasound systems connected to remote customer services receive security updates when they become available.

Recommended Actions for Your Hospital, Clinic, or Facility

Recommendations	Location of Additional Information
▪ Ensure all removable storage media is free of malware. ▪ Report any malware intrusion of the ultrasound system to the network administrator at your hospital, clinic, or facility and your service representative. ▪ Schedule periodic maintenance with your service representative.	Refer to the policies of your hospital, clinic, or facility.

Hardening the Ultrasound System

Basic protection ensures the ultrasound system is not compromised while maintaining the intended use of the ultrasound system.

- Each ultrasound system includes a BIOS/UEFI (Basic Input-Output System/Unified Extended Firmware Interface) password.
 - Change the default password after system installation for increased product security. The user is responsible to manage the changed password.
- A software-based firewall is configured with restrictions and activated.
- All user accounts not required for the intended use of the ultrasound system are disabled or deleted.
- Shared resources not required for the intended use of the medical device are disabled.
- Unnecessary network ports are closed (disabled).
- Unnecessary services are disabled.
- Unnecessary software applications are not included on the ultrasound system.
- Unauthorized software or hardware cannot be installed on the ultrasound system.
- Components on the ultrasound system are sustainable within their life cycle according to internal quality systems and external regulations for medical device manufacturers.
- The ultrasound system provides anti-virus protection by using preconfigured application allowlisting.

See also: Anti-Virus Protection, page C-2

Recommended Actions for Your Hospital, Clinic, or Facility

Recommendations	Location of Additional Information
▪ Install the ultrasound system in a secure physical location. ▪ Install the ultrasound system behind a firewall in your hospital, clinic, or facility. ▪ Store removable media in a secure location.	Refer to the policies of your hospital, clinic, or facility.
(Available only for the ultrasound system administrator) Configure the firewall settings and data export capabilities, including DICOM devices and shared network locations.	Chapter 1, System Reference Chapter 2, System Reference
(Available only for the ultrasound system administrator) Enable the Federal Information Processing Standard (FIPS) 140-2 in the configuration settings for the ultrasound system.	Chapter 6, System Reference
Configure rigorous password requirements.	Chapter 6, System Reference

Managing User Accounts and Configuring Roles

The security package on the ultrasound system provides configuration settings for managing user accounts and establishing data access permissions and functional privileges. Users have access to patient data and functions only as necessary to perform the tasks required by your hospital, clinic, or facility.

- A user login is required to access the ultrasound system. You must use the ultrasound system with the login function enabled to provide protection from security threats and vulnerabilities.
- The ultrasound system authenticates each user with a unique user name and password.
- Only personnel assigned to an administrator account can define (add, edit, or delete) user groups and user roles. The groups and roles limit access to patient data.
- Administrator privileges are not required and should not be used to operate the ultrasound system.

Recommended Actions for Your Hospital, Clinic, or Facility

Recommendations	Location of Additional Information
Assign an authorized user to an administrator account for configuring and modifying security features on the ultrasound system.	Chapter 6, System Reference
Assign the group and role to restrict access to functions and patient data on the ultrasound system.	Chapter 6, System Reference
(Available only for the ultrasound system administrator) Configure user accounts to restrict user access to the ultrasound system by requiring unique user names and passwords.	Chapter 6, System Reference
Configure rigorous password requirements. ▪ Require complex passwords. ▪ Establish periodic expirations of the passwords. ▪ Establish the quantity of unsuccessful login attempts prior to locking out a user.	Chapter 6, System Reference Refer to the policies of your hospital, clinic, or facility.
Maintain user accounts for only those who need access to the ultrasound system. Remove authorization for users who no longer require access to the ultrasound system.	Refer to the policies of your hospital, clinic, or facility.
Ensure users do not share user names and passwords.	Refer to the policies of your hospital, clinic, or facility.

Accessing the Ultrasound System for Emergency Use

Users can perform an emergency exam on the ultrasound system with an emergency login account. Access to all aspects of the system other than the functions required to perform the exam are restricted.

Recommended Actions for Your Hospital, Clinic, or Facility

Recommendations	Location of Additional Information
(Available only for the ultrasound system administrator) - Define and configure an emergency account and define access rights for an emergency user. - Track the use of the emergency user account with the audit trail logs.	Chapter 6, System Reference
- Maintain control of the emergency user account and monitor usage of the account. - Inform authorized users about the policies for utilizing an emergency user account.	Refer to the policies of your hospital, clinic, or facility.

Preventing Physical Access

The ultrasound system includes an automatic log off feature to prevent physical access by unauthorized users when the system is left idle for a period of time.

You also cannot remove device components, other than removable media, from the ultrasound system without specific tools. This design prevents physical access by unauthorized users.

Your hospital, clinic, or facility can provide additional safeguards against physical access by unauthorized users.

Recommended Actions for Your Hospital, Clinic, or Facility

Recommendations	Location of Additional Information
- Install the ultrasound system in a secure location. - Store removable media in a secure location.	Refer to the policies of your hospital, clinic, or facility.
Position the ultrasound system to prevent the viewing of patient data by unauthorized users.	Chapter 3, Instructions for Use
Configure the length of time the ultrasound system is inactive before activating the screen saver.	Chapter 1, System Reference
- Manually lock the screen to restrict physical access to the ultrasound system. A locked screen requires a user name and password to access the system. - Train authorized personnel in your hospital, clinic, or facility to use the lock screen feature when the ultrasound system is powered on and unattended.	Chapter 4, Instructions for Use

Removing or Hiding Protected Health Information

The ultrasound system is designed with features to hide protected health information.

Removable media containing patient data is under the control of your hospital, clinic, or facility.

Recommended Actions for Your Hospital, Clinic, or Facility

Recommendations	Location of Additional Information
Temporarily hide protected health information on the screen during a study.	Chapter 4, Instructions for Use
Hide protected health information from the patient report	Chapter 3, System Reference
Remove (anonymize) protected health information from a study, as required, before exporting patient data.	Chapter 5, Instructions for Use Appendix B, Instructions for Use Chapter 5, System Reference

Assuring Integrity and Confidentiality during Data Transmission

The ultrasound system is designed to ensure the integrity and confidentiality of transmitted patient data.

- Supports secure wireless network communication with the ultrasound system by authorizing a client certificate.
- Removable media containing patient data is under the control of your hospital, clinic, or facility.

Recommended Actions for Your Hospital, Clinic, or Facility

Recommendations	Location of Additional Information
Send and receive patient data using a secure wireless network transmission.	Chapter 1, System Reference
▪ Ensure the network used to transfer patient data is not accessible by non-authorized persons. ▪ Do not send patient data over an insecure network or to an insecure network device.	Refer to the policies of your hospital, clinic, or facility.

Backing Up Data and Disaster Recovery

The ultrasound system provides capabilities to recover system data and patient data. The ultrasound system is not intended as a storage device for patient data.

To prevent the loss of data, routinely back up data and settings from the ultrasound system to removable media.

Recommended Actions for Your Hospital, Clinic, or Facility

Recommendations	Location of Additional Information
Minimize the amount of patient data stored on the ultrasound system. ▪ Archive patient data to a remote DICOM archive or PACS device following acquisition of the data. ▪ Back up locally saved patient data to a DICOM archive or PACS device.	Chapter 2, System Reference Chapter 5, System Reference
Back up and store the settings of the ultrasound system on external media. Your ultrasound system administrator can subsequently restore the saved system settings.	Backing up and restoring settings is available in the service section of the configuration settings on the ultrasound system.
(Available only for the ultrasound system administrator) ▪ Archive audit logs to an external device or server. ▪ Back up and delete audit logs on a routine basis.	Chapter 6, System Reference
Remove all patient data prior to decommissioning a medical device.	Refer to the policies of your hospital, clinic, or facility.

Configuring Audit Controls and Managing Log Files

The security features of the ultrasound system provide an audit of activity on the ultrasound system by user, date, and time. The ultrasound system generates an audit log and encrypts the information in the audit log.

Recommended Actions for Your Hospital, Clinic, or Facility

Recommendations	Location of Additional Information
(Available only for the ultrasound system administrator) Configure audit controls and specify the user activities to track in an audit log.	Chapter 6, System Reference
(Available only for the ultrasound system administrator) Manage the log files on the ultrasound system. Archive and delete the log files on a regular basis.	Chapter 6, System Reference
(Available only for the ultrasound system administrator) Archive the audit logs to a removable storage device or to a device on the network; for example, a server.	Chapter 6, System Reference
Routinely review the audit logs according to the policies of your hospital, clinic, or facility.	Refer to the policies of your hospital, clinic, or facility.

Third-Party Components

The ultrasound system is maintainable and supported during the estimated life of the product.

See also: For information about open source software included on the ultrasound system, refer to the security checklist in Chapter 1 in the System Reference.

See also: For information about third-party applications, refer to the security whitepaper available from your Siemens Healthineers representative.

Remote Services

Remote customer services offers product updates of relevant operating system and product security patches.

Remote customer services encrypts all communication to and from the ultrasound system and provides the following services for your system through a virtual private network.

- Monitoring the ultrasound system by periodically collecting system properties
- Retrieving log files
- Enabling remote sessions for troubleshooting
- Auditing and documenting all changes to critical system files

Recommended Actions for Your Hospital, Clinic, or Facility

Recommendations	Location of Additional Information
▪ Immediately contact your service representative if you detect any vulnerability with the ultrasound system. ▪ Accept remote product updates.	Contact your service representative for additional information about the remote connectivity capability.

Appendix D Acoustic Output Reference

Display Resolution and Measurement Accuracy	2
Track 3, FDA 510(k) and IEC 60601-2-37 Acoustic Output Reporting	3
Summary Table for Acoustic Output	3
Definitions	5
Legend	6

Display Resolution and Measurement Accuracy

When a transducer is capable of exceeding a mechanical index (MI) of 0.01 or thermal index (TI) of 0.1, the ultrasound imaging system displays indices starting from 0.01 (MI) or 0.1 (TI) in increments of 0.01 (MI) or 0.1 (TI) for all displayed values.

A number of factors contribute to the estimation error for the displayed index. Variation among transducers and systems, approximation for real-time index calculations, and measurement errors contribute to the index display error. The measurement uncertainty for acoustic pressure, power, and center frequency is within 12%, 12%, and 7% respectively. The total estimated display accuracy is +/-15% for MI and +/-30% for TI. Definitions for these parameters can be found in the AIUM/NEMA document entitled *Standard for Real-Time Display of Thermal and Mechanical Acoustic Output Indices on Diagnostic Ultrasound Equipment* (also known as the Output Display Standard).

Track 3, FDA 510(k) and IEC 60601-2-37 Acoustic Output Reporting

Data presented in Track 3 and IEC 60601-2-37 format represents the average MI/TI values of each transducer model measured under worst-case acoustic output conditions. The on-screen MI/TI values are based on measurements on one transducer rounded up to the nearest display increment. It is possible that the values displayed on screen may exceed the MI/TI values presented in the Track 3 and IEC 60601-2-37 format.

Summary Table for Acoustic Output

An "X" indicates that either the MI index or TI indices is greater than 1.0 for each transducer/mode. A Track 3 and IEC 60601-2-37 format acoustic output table is supplied for each transducer/mode combination marked with an "X."

Transducer Model	Operating Mode							
	2D-mode	M-mode	Pulsed Wave Doppler	Color Flow Doppler	Steerable Continuous Wave Doppler	Continuous Wave Doppler	Elasticity	Contrast
(Available on the ACUSON Juniper system and the ACUSON Juniper Select system)								
5C1	X	X	X	X			X	
12L3	X	X	X	X				
11L4	X	X	X	X				
L10-5v	X	X	X	X				
5P1	X	X	X	X	X			
P4-2	X	X	X	X	X			
10MC3	X	X	X	X				
9VC2	X	X	X	X				
8VC3	X	X	X	X				
CW2						X		
CW5						X		
(Available only on the ACUSON Juniper system)								
6C1	X	X	X	X			X	X
7C2	X	X	X	X				
13L4	X	X	X	X				
14L4	X	X	X	X				X
16L4	X	X	X	X				
18H5	X	X	X	X				
11M3	X	X	X	X				
8V4	X	X	X	X	X			
10V4	X	X	X	X	X			
5VT	X	X	X	X				
9MC3	X	X	X	X				
10BCC3		X	X	X				
9VE3	X	X	X	X				
9VE4	X	X	X	X				
CW8						X		

The following rules apply to the summary table:

2D-mode	No other mode active.
M-mode	May include simultaneous 2D-mode.
Pulsed Wave Doppler Steerable Continuous Wave Doppler Continuous Wave Doppler	In duplex modes, the largest displayed TIS (scanned or non-scanned) is reported if it is larger than 1.0.
Color Flow Doppler	Includes simultaneous color flow, M-mode, 2D-mode, and Pulsed Wave Doppler. In combined modes, the largest displayed TIS (scanned or non-scanned) is reported if it is larger than 1.0.
Combined modes	The output is reported as a separate mode if the largest formulation of TIS, TIB, or TIC (if an intended use) is greater than the corresponding value reported for all constituent modes.
Other	The output is reported as a separate mode if the largest formulation of TIS, TIB, or TIC (if an intended use) is greater than the corresponding value reported for all constituent modes. TIC is reported if the transducer is intended for transcranial or neonatal cephalic use.

Definitions

Symbol		Definition		Units
FDA	IEC	FDA	IEC	
MI	<i>MI</i>	Mechanical Index	Same as FDA	n/a
TIS	<i>TIS</i>	Soft Tissue Thermal Index	Same as FDA	n/a
TIB	<i>TIB</i>	Bone Thermal Index	Same as FDA	n/a
TIC	<i>TIC</i>	Cranial Thermal Index	Same as FDA	n/a
-	$p_{ra,\alpha}$ at z_{MI}	-	Attenuated peak-rarefactional acoustic pressure at depth for mechanical index	MPa
W_o	P	Ultrasonic power, except for TIS_{scan} in which case it is the ultrasonic power passing through a one centimeter window	Output power	mW
-	P_{1x1}	-	Bounded output power	mW
z_1	z_s	Axial distance corresponding to the location of the max [min ($W_{.3}(z)$, $I_{TA,.3}(z) \times 1 \text{ cm}^2$)], where $z > z_{bp}$	Depth for TIS	cm
z_{sp}	z_b	Axial distance at which TIB is a maximum (i.e., $z_{sp} = z_{B,.3}$)	Depth for TIB	cm
-	z_{MI}	-	Depth for mechanical index	cm
$z@P_{II,.3max}$	$z_{pii,\alpha}$	The axial distance corresponding to the maximum of the derated spatial-peak pulse intensity integral (megapascals)	Depth at the point where the free-field, attenuated pulse intensity integral is a maximum	cm
f_c	f_{awf}	Center frequency	Acoustic working frequency	MHz
PRF	prf	Pulse repetition frequency	Pulse repetition rate	Hz
-	srr	-	Scan repetition rate	Hz
-	n_{pps}	-	Number of pulses per ultrasonic scan line	-
$I_{PA,.3}@P_{II,.3max}$	$I_{pa,\alpha}$ at $z_{pii,\alpha}$	Derated pulse-average intensity at the maximum of the derated spatial-peak pulse intensity integral	Same as FDA	W/cm ²
-	$I_{spta,\alpha}$ at $z_{pii,\alpha}$ OR $z_{sii,\alpha}$	-	Attenuated spatial-peak temporal-average intensity at depth for peak attenuated pulse intensity integral	mW/cm ²
-	I_{spta} at z_{pii} OR z_{sii}	-	Spatial-peak temporal-average intensity at depth for peak pulse intensity integral	mW/cm ²
$p_r@P_{II,max}$	p_r at z_{pii}	Peak rarefactional pressure at the point where the free-field, spatial-peak pulse intensity integral is a maximum	Same as FDA	MPa

Note: Pulse-average intensity (I_{pa}) is the ratio of the pulse-intensity integral (I_{pi}) to the pulse duration (t_d).

Legend

English	Translations for languages other than English
Acoustic Output Reporting Table – Track 3, FDA 510(k) and IEC 60601-2-37	n/a
(Per transducer/mode that exceeds MI or TI value of 1.0)	n/a
Transducer Model	n/a
Operating Mode	n/a
Acoustic Parameters	n/a
Index label	n/a
At surface	n/a
Below surface	n/a
units	n/a
Maximum Index Value	n/a
Index Component Value	n/a
Other Information	n/a
Operating Control Conditions	n/a
TX Level (transmit level)	n/a
Focus	n/a
a This Index is not relevant to this operating mode.	n/a
b This transducer is not intended for transcranial or neonatal cephalic uses.	n/a
c This formulation for TIS is less than that for an alternate formulation in this mode.	n/a
# No data is provided for this operation condition since the maximum index value is not reported for the reason listed.	n/a

Acoustic Output Reporting Table — Track 3, FDA 510(k) and IEC 60601-2-37

(Per transducer/mode that exceeds MI or TI value of 1.0)

Transducer model: **5C1**Operating mode: **2D-mode**

Index Label			MI	TIS		TIB		TIC
				At surface	Below surface	At surface	Below surface	
Maximum Index Value			1.50	1.71		1.71		5.40
Index Component Value				1.71	1.71	1.71	1.71	
Acoustic Parameters	IEC	Units						
	$p_{ra,\alpha}$ at z_{MI}	(MPa)	2.18					
	P	(mW)		#		187.16		712.35
	$P_{1\times1}$	(mW)		106.64		106.64		
	z_s	(cm)			#			
	z_b	(cm)					#	
	z_{MI}	(cm)	4.10					
	$z_{pii,\alpha}$	(cm)	4.10					
	f_{awf}	(MHz)	2.19	3.37		3.37		2.53
Other Information	pr_r	(Hz)	#					
	sr_r	(Hz)	28.53					
	n_{pps}	-	2.00					
	$I_{pa,\alpha}$ at $z_{pii,\alpha}$	(W/cm ²)	148.94					
	$I_{spta,\alpha}$ at $z_{pii,\alpha}$ or $z_{sii,\alpha}$	(mW/cm ²)	22.56					
	I_{spta} at z_{pii} or z_{sii}	(mW/cm ²)	41.96					
	p_r at z_{pii}	(MPa)	2.94					
Operating Control Conditions	TX Level	(dB)	0.00	0.00		0.00		0.00
	Focus	(mm)	120.00	40.00		40.00		200.00
a	This Index is not relevant to this operating mode.							
b	This transducer is not intended for transcranial or neonatal cephalic uses.							
c	This formulation for TIS is less than that for an alternate formulation in this mode.							
#	No data is provided for this operation condition since the maximum index value is not reported for the reason listed.							

Acoustic Output Reporting Table — Track 3, FDA 510(k) and IEC 60601-2-37

(Per transducer/mode that exceeds MI or TI value of 1.0)

Transducer model: **5C1**Operating mode: **M-mode**

Index Label			MI	TIS		TIB		TIC
				At surface	Below surface	At surface	Below surface	
Maximum Index Value			1.50	1.63		3.77		4.09
Index Component Value				1.16	1.63	0.83	3.77	
Acoustic Parameters	IEC	Units						
	$p_{ra,\alpha}$ at z_{MI}	(MPa)	2.18					
	P	(mW)		#		195.11		450.71
	P_{1x1}	(mW)		99.18		#		
	z_s	(cm)			4.30			
	z_b	(cm)					4.30	
	z_{MI}	(cm)	4.10					
	$z_{pii,\alpha}$	(cm)	4.10					
	f_{awf}	(MHz)	2.19	2.46		1.79		2.53
Other Information	p_{rr}	(Hz)	720.00					
	s_{rr}	(Hz)	#					
	n_{pps}	-	#					
	$I_{pa,\alpha}$ at $z_{pii,\alpha}$	(W/cm ²)	148.94					
	$I_{spta,\alpha}$ at $z_{pii,\alpha}$ or $z_{sii,\alpha}$	(mW/cm ²)	101.12					
	I_{spta} at z_{pii} or z_{sii}	(mW/cm ²)	192.78					
	p_r at z_{pii}	(MPa)	2.94					
Operating Control Conditions	TX Level	(dB)	0.00	0.00		0.00		0.00
	Focus	(mm)	120.00	100.00		60.00		200.00
a This Index is not relevant to this operating mode.								
b This transducer is not intended for transcranial or neonatal cephalic uses.								
c This formulation for TIS is less than that for an alternate formulation in this mode.								
# No data is provided for this operation condition since the maximum index value is not reported for the reason listed.								

Acoustic Output Reporting Table — Track 3, FDA 510(k) and IEC 60601-2-37

(Per transducer/mode that exceeds MI or TI value of 1.0)

Transducer model: **5C1**Operating mode: **Pulsed Wave Doppler**

Index Label			MI	TIS		TIB		TIC
				At surface	Below surface	At surface	Below surface	
Maximum Index Value			1.46	1.23		3.52		2.89
Index Component Value				0.51	1.23	0.97	3.52	
Acoustic Parameters	IEC	Units						
	$p_{ra,\alpha}$ at z_{MI}	(MPa)	2.35					
	P	(mW)		#		141.54		269.61
	P_{1x1}	(mW)		52.94		#		
	z_s	(cm)			4.60			
	z_b	(cm)					3.90	
	z_{MI}	(cm)	3.60					
	$z_{pII,\alpha}$	(cm)	3.60					
	f_{awf}	(MHz)	2.59	2.00		2.01		2.00
Other Information	p_{rr}	(Hz)	610.35					
	s_{rr}	(Hz)	#					
	n_{pps}	-	#					
	$I_{pa,\alpha}$ at $z_{pII,\alpha}$	(W/cm ²)	228.42					
	$I_{spta,\alpha}$ at $z_{pII,\alpha}$ or $z_{sII,\alpha}$	(mW/cm ²)	206.27					
	I_{spta} at z_{pII} or z_{sII}	(mW/cm ²)	386.72					
	p_r at z_{pII}	(MPa)	2.94					
Operating Control Conditions	TX Level	(dB)	0.00	0.00		0.00		0.00
	Focus	(mm)	40.00	150.00		40.00		150.00
a	This Index is not relevant to this operating mode.							
b	This transducer is not intended for transcranial or neonatal cephalic uses.							
c	This formulation for TIS is less than that for an alternate formulation in this mode.							
#	No data is provided for this operation condition since the maximum index value is not reported for the reason listed.							

Acoustic Output Reporting Table — Track 3, FDA 510(k) and IEC 60601-2-37

(Per transducer/mode that exceeds MI or TI value of 1.0)

Transducer model: **5C1**Operating mode: **Color Flow Doppler**

Index Label			MI	TIS		TIB		TIC
				At surface	Below surface	At surface	Below surface	
Maximum Index Value			1.50	1.29		1.29		4.34
Index Component Value				0.65	0.65	0.65	0.65	
Acoustic Parameters	IEC	Units						
	$p_{ra,\alpha}$ at z_{MI}	(MPa)	2.22					
	P	(mW)		#		100.56		221.39
	$P_{1\times1}$	(mW)		52.44		52.44		
	z_s	(cm)			#			
	z_b	(cm)					#	
	z_{MI}	(cm)	6.10					
	$z_{pII,\alpha}$	(cm)	6.10					
	f_{awf}	(MHz)	2.25	2.59		2.59		2.53
Other Information	p_{rr}	(Hz)	#					
	s_{rr}	(Hz)	5.89					
	n_{pps}	-	1.00					
	$I_{pa,\alpha}$ at $z_{pII,\alpha}$	(W/cm ²)	169.22					
	$I_{spta,\alpha}$ at $z_{pII,\alpha}$ or $z_{sII,\alpha}$	(mW/cm ²)	3.95					
	I_{spta} at z_{pII} or z_{sII}	(mW/cm ²)	7.45					
	p_r at z_{pII}	(MPa)	3.33					
Operating Control Conditions	TX Level	(dB)	0.00	0.00		0.00		0.00
	Focus	(mm)	80.00	40.00		40.00		120.00
a	This Index is not relevant to this operating mode.							
b	This transducer is not intended for transcranial or neonatal cephalic uses.							
c	This formulation for TIS is less than that for an alternate formulation in this mode.							
#	No data is provided for this operation condition since the maximum index value is not reported for the reason listed.							

Acoustic Output Reporting Table — Track 3, FDA 510(k) and IEC 60601-2-37

(Per transducer/mode that exceeds MI or TI value of 1.0)

Transducer model: **5C1**Operating mode: **Elasticity**

Index Label			MI	TIS		TIB		TIC
				At surface	Below surface	At surface	Below surface	
Maximum Index Value			1.49	0.90		2.21		1.97
Index Component Value				0.36	0.78	0.39	2.08	
Acoustic Parameters	IEC	Units						
	$p_{ra,\alpha}$ at z_{MI}	(MPa)	2.58					
	P	(mW)		#		128.27		141.96
	$P_{1\times1}$	(mW)		37.64		#		
	z_s	(cm)			4.00			
	z_b	(cm)					6.20	
	z_{MI}	(cm)	3.50					
	$z_{pii,\alpha}$	(cm)	3.50					
	f_{awf}	(MHz)	3.05	2.02		2.03		2.02
Other Information	p_{rr}	(Hz)	39.10					
	s_{rr}	(Hz)	#					
	n_{pps}	-	#					
	$I_{pa,\alpha}$ at $z_{pii,\alpha}$	(W/cm ²)	359.31					
	$I_{spta,\alpha}$ at $z_{pii,\alpha}$ or $z_{sii,\alpha}$	(mW/cm ²)	27.42					
	I_{spta} at z_{pii} or z_{sii}	(mW/cm ²)	43.54					
	p_r at z_{pii}	(MPa)	3.42					
Operating Control Conditions	TX Level	(dB)	0.00	0.00		0.00		0.00
	Focus	(mm)	50.00	92.00		77.00		92.00
a	This Index is not relevant to this operating mode.							
b	This transducer is not intended for transcranial or neonatal cephalic uses.							
c	This formulation for TIS is less than that for an alternate formulation in this mode.							
#	No data is provided for this operation condition since the maximum index value is not reported for the reason listed.							

Acoustic Output Reporting Table — Track 3, FDA 510(k) and IEC 60601-2-37

(Per transducer/mode that exceeds MI or TI value of 1.0)

Transducer model: **6C1**Operating mode: **2D-mode**

Index Label			MI	TIS		TIB		TIC
				At surface	Below surface	At surface	Below surface	
Maximum Index Value			1.39	1.87		1.87		5.31
Index Component Value				1.87	1.87	1.87	1.87	
Acoustic Parameters	IEC	Units						
	$p_{ra,\alpha}$ at z_{MI}	(MPa)	2.04					
	P	(mW)		#		167.79		626.88
	P_{1x1}	(mW)		112.35		112.35		
	z_s	(cm)			#			
	z_b	(cm)					#	
	z_{MI}	(cm)	4.00					
	$z_{pii,\alpha}$	(cm)	4.00					
	f_{awf}	(MHz)	2.19	3.49		3.49		1.76
Other Information	p_{rr}	(Hz)	#					
	s_{rr}	(Hz)	9.30					
	n_{pps}	-	2.00					
	$I_{pa,\alpha}$ at $z_{pii,\alpha}$	(W/cm ²)	322.37					
	$I_{spta,\alpha}$ at $z_{pii,\alpha}$ or $z_{sii,\alpha}$	(mW/cm ²)	4.69					
	I_{spta} at z_{pii} or z_{sii}	(mW/cm ²)	8.71					
	p_r at z_{pii}	(MPa)	2.69					
Operating Control Conditions	TX Level	(dB)	0.00	0.00		0.00		0.00
	Focus	(mm)	60.00	40.00		40.00		190.00
a This Index is not relevant to this operating mode.								
b This transducer is not intended for transcranial or neonatal cephalic uses.								
c This formulation for TIS is less than that for an alternate formulation in this mode.								
# No data is provided for this operation condition since the maximum index value is not reported for the reason listed.								

Acoustic Output Reporting Table — Track 3, FDA 510(k) and IEC 60601-2-37

(Per transducer/mode that exceeds MI or TI value of 1.0)

Transducer model: **6C1**Operating mode: **M-mode**

Index Label			MI	TIS		TIB		TIC
				At surface	Below surface	At surface	Below surface	
Maximum Index Value			1.39	1.91		5.36		4.92
Index Component Value				0.96	1.91	0.96	5.36	
Acoustic Parameters	IEC	Units						
	$p_{ra,\alpha}$ at z_{MI}	(MPa)	1.97					
	P	(mW)		#		377.38		377.38
	$P_{1\times1}$	(mW)		130.58		#		
	z_s	(cm)			3.50			
	z_b	(cm)					3.80	
	z_{MI}	(cm)	3.80					
	$z_{pII,\alpha}$	(cm)	3.80					
	f_{awf}	(MHz)	2.07	1.54		1.54		1.54
Other Information	prr	(Hz)	1440.09					
	srr	(Hz)	#					
	n_{pps}	-	#					
	$I_{pa,\alpha}$ at $z_{pII,\alpha}$	(W/cm ²)	197.66					
	$I_{spta,\alpha}$ at $z_{pII,\alpha}$ or $z_{sII,\alpha}$	(mW/cm ²)	192.02					
	I_{spta} at z_{pII} or z_{sII}	(mW/cm ²)	337.61					
	p_r at z_{pII}	(MPa)	2.36					
	Operating Control Conditions	TX Level	(dB)	0.00	0.00		0.00	
Focus		(mm)	240.00	190.00		190.00		190.00
a	This Index is not relevant to this operating mode.							
b	This transducer is not intended for transcranial or neonatal cephalic uses.							
c	This formulation for TIS is less than that for an alternate formulation in this mode.							
#	No data is provided for this operation condition since the maximum index value is not reported for the reason listed.							

Acoustic Output Reporting Table — Track 3, FDA 510(k) and IEC 60601-2-37

(Per transducer/mode that exceeds MI or TI value of 1.0)

Transducer model: **6C1**Operating mode: **Pulsed Wave Doppler**

Index Label			MI	TIS		TIB		TIC
				At surface	Below surface	At surface	Below surface	
Maximum Index Value			1.39	2.43		5.33		4.47
Index Component Value				1.23	2.43	1.07	5.33	
Acoustic Parameters	IEC	Units						
	$p_{ra,\alpha}$ at z_{MI}	(MPa)	2.25					
	P	(mW)		#		358.45		450.16
	P_{1x1}	(mW)		99.05		#		
	z_s	(cm)			4.00			
	z_b	(cm)					4.30	
	z_{MI}	(cm)	5.20					
	$z_{pii,\alpha}$	(cm)	5.20					
	f_{awf}	(MHz)	2.64	2.61		2.05		2.04
Other Information	p_{rr}	(Hz)	610.35					
	s_{rr}	(Hz)	#					
	n_{pps}	-	#					
	$I_{pa,\alpha}$ at $z_{pii,\alpha}$	(W/cm ²)	341.99					
	$I_{spta,\alpha}$ at $z_{pii,\alpha}$ or $z_{sii,\alpha}$	(mW/cm ²)	304.24					
	I_{spta} at z_{pii} or z_{sii}	(mW/cm ²)	788.66					
	p_r at z_{pii}	(MPa)	3.37					
Operating Control Conditions	TX Level	(dB)	0.00	0.00		0.00		0.00
	Focus	(mm)	60.00	100.00		80.00		120.00
a This Index is not relevant to this operating mode.								
b This transducer is not intended for transcranial or neonatal cephalic uses.								
c This formulation for TIS is less than that for an alternate formulation in this mode.								
# No data is provided for this operation condition since the maximum index value is not reported for the reason listed.								

Acoustic Output Reporting Table — Track 3, FDA 510(k) and IEC 60601-2-37

(Per transducer/mode that exceeds MI or TI value of 1.0)

Transducer model: **6C1**Operating mode: **Color Flow Doppler**

Index Label			MI	TIS		TIB		TIC
				At surface	Below surface	At surface	Below surface	
Maximum Index Value			1.39	2.23		2.23		6.00
Index Component Value				2.14	2.14	2.14	2.14	
Acoustic Parameters	IEC	Units						
	$p_{ra,\alpha}$ at z_{MI}	(MPa)	1.99					
	P	(mW)		#		243.10		711.19
	$P_{1\times1}$	(mW)		135.59		135.59		
	z_s	(cm)			#			
	z_b	(cm)					#	
	z_{MI}	(cm)	4.90					
	$z_{pii,\alpha}$	(cm)	4.90					
	f_{awf}	(MHz)	2.05	3.32		3.32		2.05
Other Information	p_{rr}	(Hz)	#					
	s_{rr}	(Hz)	140.93					
	n_{pps}	-	9.00					
	$I_{pa,\alpha}$ at $z_{pii,\alpha}$	(W/cm ²)	361.65					
	$I_{spta,\alpha}$ at $z_{pii,\alpha}$ or $z_{sii,\alpha}$	(mW/cm ²)	161.11					
	I_{spta} at z_{pii} or z_{sii}	(mW/cm ²)	275.89					
	p_r at z_{pii}	(MPa)	2.74					
Operating Control Conditions	TX Level	(dB)	0.00	0.00		0.00		0.00
	Focus	(mm)	60.00	40.00		40.00		240.00
a	This Index is not relevant to this operating mode.							
b	This transducer is not intended for transcranial or neonatal cephalic uses.							
c	This formulation for TIS is less than that for an alternate formulation in this mode.							
#	No data is provided for this operation condition since the maximum index value is not reported for the reason listed.							

Acoustic Output Reporting Table — Track 3, FDA 510(k) and IEC 60601-2-37

(Per transducer/mode that exceeds MI or TI value of 1.0)

Transducer model: **6C1**Operating mode: **Elasticity**

Index Label			MI	TIS		TIB		TIC
				At surface	Below surface	At surface	Below surface	
Maximum Index Value			1.51	1.44		3.10		3.40
Index Component Value				0.81	1.27	0.82	2.92	
Acoustic Parameters	IEC	Units						
	$p_{ra,\alpha}$ at z_{MI}	(MPa)	2.39					
	P	(mW)		#		204.72		289.24
	$P_{1\times1}$	(mW)		67.96		#		
	z_s	(cm)			4.20			
	z_b	(cm)					5.60	
	z_{MI}	(cm)	3.90					
	$z_{pII,\alpha}$	(cm)	3.90					
	f_{awf}	(MHz)	2.50	2.50		2.50		2.48
Other Information	p_{rr}	(Hz)	#					
	s_{rr}	(Hz)	2.21					
	n_{pps}	-	#					
	$I_{pa,\alpha}$ at $z_{pII,\alpha}$	(W/cm ²)	475.49					
	$I_{spta,\alpha}$ at $z_{pII,\alpha}$ or $z_{sII,\alpha}$	(mW/cm ²)	322.42					
	I_{spta} at z_{pII} or z_{sII}	(mW/cm ²)	655.26					
	p_r at z_{pII}	(MPa)	3.26					
Operating Control Conditions	TX Level	(dB)	0.00	0.00		0.00		0.00
	Focus	(mm)	42.00	66.00		60.50		92.00
a This Index is not relevant to this operating mode.								
b This transducer is not intended for transcranial or neonatal cephalic uses.								
c This formulation for TIS is less than that for an alternate formulation in this mode.								
# No data is provided for this operation condition since the maximum index value is not reported for the reason listed.								

Acoustic Output Reporting Table — Track 3, FDA 510(k) and IEC 60601-2-37

(Per transducer/mode that exceeds MI or TI value of 1.0)

Transducer model: **6C1**Operating mode: **Contrast**

Index Label			MI	TIS		TIB		TIC
				At surface	Below surface	At surface	Below surface	
Maximum Index Value			1.38	1.88		1.88		5.46
Index Component Value				1.43	1.43	1.43	1.43	
Acoustic Parameters	IEC	Units						
	$p_{ra,\alpha}$ at z_{MI}	(MPa)	1.75					
	P	(mW)		#		206.77		523.54
	$P_{1\times1}$	(mW)		122.20		122.20		
	z_s	(cm)			#			
	z_b	(cm)					#	
	z_{MI}	(cm)	6.30					
	$z_{pii,\alpha}$	(cm)	6.30					
	f_{awf}	(MHz)	1.67	2.46		2.46		2.47
Other Information	pr_r	(Hz)	#					
	srr	(Hz)	22.68					
	n_{pps}	-	2.00					
	$I_{pa,\alpha}$ at $z_{pii,\alpha}$	(W/cm ²)	159.67					
	$I_{spta,\alpha}$ at $z_{pii,\alpha}$ or $z_{sii,\alpha}$	(mW/cm ²)	46.43					
	I_{spta} at z_{pii} or z_{sii}	(mW/cm ²)	67.85					
	p_r at z_{pii}	(MPa)	2.25					
Operating Control Conditions	TX Level	(dB)	0.00	0.00		0.00		0.00
	Focus	(mm)	80.00	40.00		40.00		190.00
a	This Index is not relevant to this operating mode.							
b	This transducer is not intended for transcranial or neonatal cephalic uses.							
c	This formulation for TIS is less than that for an alternate formulation in this mode.							
#	No data is provided for this operation condition since the maximum index value is not reported for the reason listed.							

Acoustic Output Reporting Table — Track 3, FDA 510(k) and IEC 60601-2-37

(Per transducer/mode that exceeds MI or TI value of 1.0)

Transducer model: **7C2**Operating mode: **2D-mode**

Index Label			MI	TIS		TIB		TIC
				At surface	Below surface	At surface	Below surface	
Maximum Index Value			1.38	3.49		3.49		6.00
Index Component Value				3.49	3.49	3.49	3.49	
Acoustic Parameters	IEC	Units						
	$p_{ra,\alpha}$ at z_{MI}	(MPa)	2.48					
	P	(mW)		#		209.61		728.40
	P_{1x1}	(mW)		163.91		163.91		
	z_s	(cm)			#			
	z_b	(cm)					#	
	z_{MI}	(cm)	0.50					
	$z_{pii,\alpha}$	(cm)	1.30					
	f_{awf}	(MHz)	3.20	4.46		4.46		3.10
Other Information	p_{rr}	(Hz)	#					
	s_{rr}	(Hz)	65.62					
	n_{pps}	-	1.00					
	$I_{pa,\alpha}$ at $z_{pii,\alpha}$	(W/cm ²)	252.60					
	$I_{spta,\alpha}$ at $z_{pii,\alpha}$ or $z_{sii,\alpha}$	(mW/cm ²)	8.27					
	I_{spta} at z_{pii} or z_{sii}	(mW/cm ²)	20.89					
	p_r at z_{pii}	(MPa)	2.14					
Operating Control Conditions	TX Level	(dB)	0.00	0.00		0.00		0.00
	Focus	(mm)	20.00	20.00		20.00		200.00
a This Index is not relevant to this operating mode.								
b This transducer is not intended for transcranial or neonatal cephalic uses.								
c This formulation for TIS is less than that for an alternate formulation in this mode.								
# No data is provided for this operation condition since the maximum index value is not reported for the reason listed.								

Acoustic Output Reporting Table — Track 3, FDA 510(k) and IEC 60601-2-37

(Per transducer/mode that exceeds MI or TI value of 1.0)

Transducer model: **7C2**Operating mode: **M-mode**

Index Label			MI	TIS		TIB		TIC
				At surface	Below surface	At surface	Below surface	
Maximum Index Value			1.38	2.41		2.87		4.98
Index Component Value				2.41	1.84	2.41	2.87	
Acoustic Parameters	IEC	Units						
	$p_{ra,\alpha}$ at z_{MI}	(MPa)	2.48					
	P	(mW)		#		184.64		417.09
	$P_{1\times1}$	(mW)		159.38		159.38		
	z_s	(cm)			1.90			
	z_b	(cm)					2.00	
	z_{MI}	(cm)	0.50					
	$z_{pii,\alpha}$	(cm)	1.30					
	f_{awf}	(MHz)	3.20	3.18		3.18		3.15
Other Information	prr	(Hz)	720.00					
	srr	(Hz)	#					
	n_{pps}	-	#					
	$I_{pa,\alpha}$ at $z_{pii,\alpha}$	(W/cm ²)	252.60					
	$I_{spta,\alpha}$ at $z_{pii,\alpha}$ or $z_{sii,\alpha}$	(mW/cm ²)	149.95					
	I_{spta} at z_{pii} or z_{sii}	(mW/cm ²)	170.19					
	p_r at z_{pii}	(MPa)	2.14					
Operating Control Conditions	TX Level	(dB)	0.00	0.00		0.00		0.00
	Focus	(mm)	20.00	40.00		40.00		120.00
a This Index is not relevant to this operating mode.								
b This transducer is not intended for transcranial or neonatal cephalic uses.								
c This formulation for TIS is less than that for an alternate formulation in this mode.								
# No data is provided for this operation condition since the maximum index value is not reported for the reason listed.								

Acoustic Output Reporting Table — Track 3, FDA 510(k) and IEC 60601-2-37

(Per transducer/mode that exceeds MI or TI value of 1.0)

Transducer model: **7C2**Operating mode: **Pulsed Wave Doppler**

Index Label			MI	TIS		TIB		TIC
				At surface	Below surface	At surface	Below surface	
Maximum Index Value			1.49	2.09		3.80		3.62
Index Component Value				2.09	1.66	1.74	3.80	
Acoustic Parameters	IEC	Units						
	$p_{ra,\alpha}$ at z_{MI}	(MPa)	2.71					
	P	(mW)		#		110.30		360.54
	P_{1x1}	(mW)		137.84		110.30		
	z_s	(cm)			1.90			
	z_b	(cm)					1.40	
	z_{MI}	(cm)	1.20					
	$z_{pii,\alpha}$	(cm)	1.30					
	f_{awf}	(MHz)	3.31	3.18		3.31		3.16
Other Information	p_{rr}	(Hz)	1220.70					
	s_{rr}	(Hz)	#					
	n_{pps}	-	#					
	$I_{pa,\alpha}$ at $z_{pii,\alpha}$	(W/cm ²)	366.23					
	$I_{spta,\alpha}$ at $z_{pii,\alpha}$ or $z_{sii,\alpha}$	(mW/cm ²)	473.78					
	I_{spta} at z_{pii} or z_{sii}	(mW/cm ²)	637.82					
	p_r at z_{pii}	(MPa)	3.04					
Operating Control Conditions	TX Level	(dB)	0.00	0.00		0.00		0.00
	Focus	(mm)	20.00	40.00		20.00		200.00
a This Index is not relevant to this operating mode.								
b This transducer is not intended for transcranial or neonatal cephalic uses.								
c This formulation for TIS is less than that for an alternate formulation in this mode.								
# No data is provided for this operation condition since the maximum index value is not reported for the reason listed.								

Acoustic Output Reporting Table — Track 3, FDA 510(k) and IEC 60601-2-37

(Per transducer/mode that exceeds MI or TI value of 1.0)

Transducer model: **7C2**Operating mode: **Color Flow Doppler**

Index Label			MI	TIS		TIB		TIC
				At surface	Below surface	At surface	Below surface	
Maximum Index Value			1.49	0.99		0.99		0.74
Index Component Value				0.99	0.99	0.99	0.99	
Acoustic Parameters	IEC	Units						
	$p_{ra,\alpha}$ at z_{MI}	(MPa)	2.48					
	P	(mW)		#		84.38		49.11
	P_{1x1}	(mW)		65.60		65.60		
	z_s	(cm)			#			
	z_b	(cm)					#	
	z_{MI}	(cm)	1.00					
	$z_{pii,\alpha}$	(cm)	1.50					
	f_{awf}	(MHz)	2.75	3.18		3.18		3.15
Other Information	prf	(Hz)	#					
	srr	(Hz)	114.86					
	n_{pps}	-	12.00					
	$I_{pa,\alpha}$ at $z_{pii,\alpha}$	(W/cm ²)	225.09					
	$I_{spta,\alpha}$ at $z_{pii,\alpha}$ or $z_{sii,\alpha}$	(mW/cm ²)	11.85					
	I_{spta} at z_{pii} or z_{sii}	(mW/cm ²)	25.78					
	p_r at z_{pii}	(MPa)	2.73					
Operating Control Conditions	TX Level	(dB)	0.00	0.00		0.00		0.00
	Focus	(mm)	20.00	40.00		40.00		60.00
a	This Index is not relevant to this operating mode.							
b	This transducer is not intended for transcranial or neonatal cephalic uses.							
c	This formulation for TIS is less than that for an alternate formulation in this mode.							
#	No data is provided for this operation condition since the maximum index value is not reported for the reason listed.							

Acoustic Output Reporting Table — Track 3, FDA 510(k) and IEC 60601-2-37

(Per transducer/mode that exceeds MI or TI value of 1.0)

Transducer model: **12L3**Operating mode: **2D-mode**

Index Label			MI	TIS		TIB		TIC
				At surface	Below surface	At surface	Below surface	
Maximum Index Value			1.14	0.93		0.93		1.87
Index Component Value				0.93	0.93	0.93	0.93	
Acoustic Parameters	IEC	Units						
	$p_{ra,\alpha}$ at z_{MI}	(MPa)	2.27					
	P	(mW)		#		159.52		159.52
	P_{1x1}	(mW)		31.28		31.28		
	z_s	(cm)			#			
	z_b	(cm)					#	
	z_{MI}	(cm)	2.00					
	$z_{pii,\alpha}$	(cm)	2.00					
	f_{awf}	(MHz)	4.09	6.27		6.27		6.27
Other Information	prf	(Hz)	#					
	srf	(Hz)	79.52					
	n_{pps}	-	2.00					
	$I_{pa,\alpha}$ at $z_{pii,\alpha}$	(W/cm ²)	255.40					
	$I_{spta,\alpha}$ at $z_{pii,\alpha}$ or $z_{sii,\alpha}$	(mW/cm ²)	18.13					
	I_{spta} at z_{pii} or z_{sii}	(mW/cm ²)	29.37					
	p_r at z_{pii}	(MPa)	2.88					
Operating Control Conditions	TX Level	(dB)	0.00	0.00		0.00		0.00
	Focus	(mm)	20.00	80.00		80.00		80.00
a	This Index is not relevant to this operating mode.							
b	This transducer is not intended for transcranial or neonatal cephalic uses.							
c	This formulation for TIS is less than that for an alternate formulation in this mode.							
#	No data is provided for this operation condition since the maximum index value is not reported for the reason listed.							

Acoustic Output Reporting Table — Track 3, FDA 510(k) and IEC 60601-2-37

(Per transducer/mode that exceeds MI or TI value of 1.0)

Transducer model: **12L3**Operating mode: **M-mode**

Index Label			MI	TIS		TIB		TIC
				At surface	Below surface	At surface	Below surface	
Maximum Index Value			0.83	1.13		1.13		0.97
Index Component Value				1.13	0.63	1.13	0.52	
Acoustic Parameters	IEC	Units						
	$p_{ra,\alpha}$ at z_{MI}	(MPa)	2.08					
	P	(mW)		#		50.37		50.37
	$P_{1\times1}$	(mW)		37.69		37.69		
	z_s	(cm)			2.00			
	z_b	(cm)					2.00	
	z_{MI}	(cm)	1.20					
	$z_{pii,\alpha}$	(cm)	2.00					
	f_{awf}	(MHz)	6.04	6.27		6.27		6.27
Other Information	p_{rr}	(Hz)	720.00					
	s_{rr}	(Hz)	#					
	n_{pps}	-	#					
	$I_{pa,\alpha}$ at $z_{pii,\alpha}$	(W/cm ²)	95.56					
	$I_{spta,\alpha}$ at $z_{pii,\alpha}$ or $z_{sii,\alpha}$	(mW/cm ²)	26.97					
	I_{spta} at z_{pii} or z_{sii}	(mW/cm ²)	78.26					
	p_r at z_{pii}	(MPa)	2.61					
Operating Control Conditions	TX Level	(dB)	0.00	0.00		0.00		0.00
	Focus	(mm)	60.00	80.00		80.00		80.00
a	This Index is not relevant to this operating mode.							
b	This transducer is not intended for transcranial or neonatal cephalic uses.							
c	This formulation for TIS is less than that for an alternate formulation in this mode.							
#	No data is provided for this operation condition since the maximum index value is not reported for the reason listed.							

Acoustic Output Reporting Table — Track 3, FDA 510(k) and IEC 60601-2-37

(Per transducer/mode that exceeds MI or TI value of 1.0)

Transducer model: **12L3**Operating mode: **Pulsed Wave Doppler**

Index Label			MI	TIS		TIB		TIC
				At surface	Below surface	At surface	Below surface	
Maximum Index Value			1.46	1.49		1.88		1.73
Index Component Value				1.49	0.48	0.62	1.88	
Acoustic Parameters	IEC	Units						
	$p_{ra,\alpha}$ at z_{MI}	(MPa)	2.91					
	P	(mW)		#		32.59		26.13
	P_{1x1}	(mW)		50.06		32.59		
	z_s	(cm)			1.70			
	z_b	(cm)					1.50	
	z_{MI}	(cm)	2.40					
	$z_{pii,\alpha}$	(cm)	2.40					
	f_{awf}	(MHz)	4.00	6.25		4.02		4.05
Other Information	p_{rr}	(Hz)	610.35					
	s_{rr}	(Hz)	#					
	n_{pps}	-	#					
	$I_{pa,\alpha}$ at $z_{pii,\alpha}$	(W/cm ²)	592.67					
	$I_{spta,\alpha}$ at $z_{pii,\alpha}$ or $z_{sii,\alpha}$	(mW/cm ²)	372.18					
	I_{spta} at z_{pii} or z_{sii}	(mW/cm ²)	713.83					
	p_r at z_{pii}	(MPa)	3.98					
Operating Control Conditions	TX Level	(dB)	0.00	0.00		0.00		0.00
	Focus	(mm)	30.00	40.00		15.00		5.00
a This Index is not relevant to this operating mode.								
b This transducer is not intended for transcranial or neonatal cephalic uses.								
c This formulation for TIS is less than that for an alternate formulation in this mode.								
# No data is provided for this operation condition since the maximum index value is not reported for the reason listed.								

Acoustic Output Reporting Table — Track 3, FDA 510(k) and IEC 60601-2-37

(Per transducer/mode that exceeds MI or TI value of 1.0)

Transducer model: **12L3**Operating mode: **Color Flow Doppler**

Index Label			MI	TIS		TIB		TIC
				At surface	Below surface	At surface	Below surface	
Maximum Index Value			1.46	1.12		1.12		1.64
Index Component Value				1.12	1.12	1.12	1.12	
Acoustic Parameters	IEC	Units						
	$p_{ra,\alpha}$ at z_{MI}	(MPa)	2.92					
	P	(mW)		#		49.54		140.28
	P_{1x1}	(mW)		37.06		37.06		
	z_s	(cm)			#			
	z_b	(cm)					#	
	z_{MI}	(cm)	1.90					
	$z_{pII,\alpha}$	(cm)	1.90					
	f_{awf}	(MHz)	4.02	6.35		6.35		4.00
Other Information	p_{rr}	(Hz)	#					
	s_{rr}	(Hz)	68.15					
	n_{pps}	-	12.00					
	$I_{pa,\alpha}$ at $z_{pII,\alpha}$	(W/cm ²)	558.45					
	$I_{spta,\alpha}$ at $z_{pII,\alpha}$ or $z_{sII,\alpha}$	(mW/cm ²)	24.79					
	I_{spta} at z_{pII} or z_{sII}	(mW/cm ²)	42.91					
	p_r at z_{pII}	(MPa)	3.80					
Operating Control Conditions	TX Level	(dB)	0.00	0.00		0.00		0.00
	Focus	(mm)	20.00	30.00		30.00		30.00
a	This Index is not relevant to this operating mode.							
b	This transducer is not intended for transcranial or neonatal cephalic uses.							
c	This formulation for TIS is less than that for an alternate formulation in this mode.							
#	No data is provided for this operation condition since the maximum index value is not reported for the reason listed.							

Acoustic Output Reporting Table — Track 3, FDA 510(k) and IEC 60601-2-37

(Per transducer/mode that exceeds MI or TI value of 1.0)

Transducer model: **11L4**Operating mode: **2D-mode**

Index Label			MI	TIS		TIB		TIC
				At surface	Below surface	At surface	Below surface	
Maximum Index Value			1.28	2.48		2.48		2.55
Index Component Value				2.48	2.48	2.48	2.48	
Acoustic Parameters	IEC	Units						
	$p_{ra,\alpha}$ at z_{MI}	(MPa)	3.24					
	P	(mW)		#		109.47		159.33
	$P_{1\times1}$	(mW)		79.97		79.97		
	z_s	(cm)			#			
	z_b	(cm)					#	
	z_{MI}	(cm)	1.20					
	$z_{pII,\alpha}$	(cm)	1.20					
	f_{awf}	(MHz)	6.61	6.51		6.51		6.67
Other Information	p_{rr}	(Hz)	#					
	s_{rr}	(Hz)	58.09					
	n_{pps}	-	1.00					
	$I_{pa,\alpha}$ at $z_{pII,\alpha}$	(W/cm ²)	325.78					
	$I_{spta,\alpha}$ at $z_{pII,\alpha}$ or $z_{sII,\alpha}$	(mW/cm ²)	1.70					
	I_{spta} at z_{pII} or z_{sII}	(mW/cm ²)	2.94					
	p_r at z_{pII}	(MPa)	3.73					
Operating Control Conditions	TX Level	(dB)	0.00	0.00		0.00		0.00
	Focus	(mm)	30.00	50.00		50.00		70.00
a	This Index is not relevant to this operating mode.							
b	This transducer is not intended for transcranial or neonatal cephalic uses.							
c	This formulation for TIS is less than that for an alternate formulation in this mode.							
#	No data is provided for this operation condition since the maximum index value is not reported for the reason listed.							

Acoustic Output Reporting Table — Track 3, FDA 510(k) and IEC 60601-2-37

(Per transducer/mode that exceeds MI or TI value of 1.0)

Transducer model: **11L4**Operating mode: **M-mode**

Index Label			MI	TIS		TIB		TIC
				At surface	Below surface	At surface	Below surface	
Maximum Index Value			1.50	1.23		1.23		0.98
Index Component Value				1.23	0.61	1.23	0.70	
Acoustic Parameters	IEC	Units						
	$p_{ra,\alpha}$ at z_{MI}	(MPa)	3.80					
	P	(mW)		#		41.06		41.06
	P_{1x1}	(mW)		41.06		41.06		
	z_s	(cm)			1.60			
	z_b	(cm)					1.60	
	z_{MI}	(cm)	1.20					
	$z_{pII,\alpha}$	(cm)	1.20					
	f_{awf}	(MHz)	6.59	6.27		6.27		6.27
	p_{rr}	(Hz)	720.00					
Other Information	s_{rr}	(Hz)	#					
	n_{pps}	-	#					
	$I_{pa,\alpha}$ at $z_{pII,\alpha}$	(W/cm ²)	518.39					
	$I_{spta,\alpha}$ at $z_{pII,\alpha}$ or $z_{sII,\alpha}$	(mW/cm ²)	78.52					
	I_{spta} at z_{pII} or z_{sII}	(mW/cm ²)	139.42					
	p_r at z_{pII}	(MPa)	4.84					
	TX Level	(dB)	0.00	0.00		0.00		0.00
	Focus	(mm)	15.00	70.00		70.00		70.00
Operating Control Conditions								
a This Index is not relevant to this operating mode. b This transducer is not intended for transcranial or neonatal cephalic uses. c This formulation for TIS is less than that for an alternate formulation in this mode. # No data is provided for this operation condition since the maximum index value is not reported for the reason listed.								

Acoustic Output Reporting Table — Track 3, FDA 510(k) and IEC 60601-2-37

(Per transducer/mode that exceeds MI or TI value of 1.0)

Transducer model: **11L4**Operating mode: **Pulsed Wave Doppler**

Index Label			MI	TIS		TIB		TIC
				At surface	Below surface	At surface	Below surface	
Maximum Index Value			1.49	2.31		2.35		1.96
Index Component Value				2.31	1.19	2.31	2.35	
Acoustic Parameters	IEC	Units						
	$p_{ra,\alpha}$ at z_{MI}	(MPa)	3.03					
	P	(mW)		#		81.15		81.15
	P_{1x1}	(mW)		81.15		81.15		
	z_s	(cm)			1.60			
	z_b	(cm)					1.60	
	z_{MI}	(cm)	1.50					
	$z_{pII,\alpha}$	(cm)	1.50					
	f_{awf}	(MHz)	4.07	5.99		5.99		5.99
	p_{rr}	(Hz)	976.56					
Other Information	s_{rr}	(Hz)	#					
	n_{pps}	-	#					
	$I_{pa,\alpha}$ at $z_{pII,\alpha}$	(W/cm ²)	293.81					
	$I_{spta,\alpha}$ at $z_{pII,\alpha}$ or $z_{sII,\alpha}$	(mW/cm ²)	336.41					
	I_{spta} at z_{pII} or z_{sII}	(mW/cm ²)	525.61					
	p_r at z_{pII}	(MPa)	3.73					
	TX Level	(dB)	0.00	0.00		0.00		0.00
	Focus	(mm)	20.00	50.00		50.00		50.00
Operating Control Conditions								
a This Index is not relevant to this operating mode. b This transducer is not intended for transcranial or neonatal cephalic uses. c This formulation for TIS is less than that for an alternate formulation in this mode. # No data is provided for this operation condition since the maximum index value is not reported for the reason listed.								

Acoustic Output Reporting Table — Track 3, FDA 510(k) and IEC 60601-2-37

(Per transducer/mode that exceeds MI or TI value of 1.0)

Transducer model: **11L4**Operating mode: **Color Flow Doppler**

Index Label			MI	TIS		TIB		TIC
				At surface	Below surface	At surface	Below surface	
Maximum Index Value			1.50	1.60		1.60		1.02
Index Component Value				1.60	1.60	1.60	1.60	
Acoustic Parameters	IEC	Units						
	$p_{ra,\alpha}$ at z_{MI}	(MPa)	3.70					
	P	(mW)		#		59.06		49.90
	$P_{1\times1}$	(mW)		55.09		55.09		
	z_s	(cm)			#			
	z_b	(cm)					#	
	z_{MI}	(cm)	0.60					
	$z_{pII,\alpha}$	(cm)	1.00					
	f_{awf}	(MHz)	6.04	6.08		6.08		6.08
Other Information	prr	(Hz)	#					
	srr	(Hz)	169.85					
	n_{pps}	-	12.00					
	$I_{pa,\alpha}$ at $z_{pII,\alpha}$	(W/cm ²)	559.68					
	$I_{spta,\alpha}$ at $z_{pII,\alpha}$ or $z_{sII,\alpha}$	(mW/cm ²)	93.54					
	I_{spta} at z_{pII} or z_{sII}	(mW/cm ²)	147.90					
	p_r at z_{pII}	(MPa)	4.09					
Operating Control Conditions	TX Level	(dB)	0.00	0.00		0.00		0.00
	Focus	(mm)	10.00	7.50		7.50		7.50
a	This Index is not relevant to this operating mode.							
b	This transducer is not intended for transcranial or neonatal cephalic uses.							
c	This formulation for TIS is less than that for an alternate formulation in this mode.							
#	No data is provided for this operation condition since the maximum index value is not reported for the reason listed.							

Acoustic Output Reporting Table — Track 3, FDA 510(k) and IEC 60601-2-37

(Per transducer/mode that exceeds MI or TI value of 1.0)

Transducer model: **13L4**Operating mode: **2D-mode**

Index Label			MI	TIS		TIB		TIC
				At surface	Below surface	At surface	Below surface	
Maximum Index Value			1.42	1.04		1.04		1.12
Index Component Value				1.04	1.04	1.04	1.04	
Acoustic Parameters	IEC	Units						
	$p_{ra,\alpha}$ at z_{MI}	(MPa)	3.70					
	P	(mW)		#		34.60		66.55
	P_{1x1}	(mW)		32.41		32.41		
	z_s	(cm)			#			
	z_b	(cm)					#	
	z_{MI}	(cm)	0.60					
	$z_{pii,\alpha}$	(cm)	0.70					
	f_{awf}	(MHz)	6.73	6.76		6.76		6.36
Other Information	p_{rr}	(Hz)	#					
	s_{rr}	(Hz)	98.89					
	n_{pps}	-	2.00					
	$I_{pa,\alpha}$ at $z_{pii,\alpha}$	(W/cm ²)	473.70					
	$I_{spta,\alpha}$ at $z_{pii,\alpha}$ or $z_{sii,\alpha}$	(mW/cm ²)	12.64					
	I_{spta} at z_{pii} or z_{sii}	(mW/cm ²)	18.32					
	p_r at z_{pii}	(MPa)	4.27					
Operating Control Conditions	TX Level	(dB)	0.00	0.00		0.00		0.00
	Focus	(mm)	10.00	5.00		5.00		50.00
a This Index is not relevant to this operating mode.								
b This transducer is not intended for transcranial or neonatal cephalic uses.								
c This formulation for TIS is less than that for an alternate formulation in this mode.								
# No data is provided for this operation condition since the maximum index value is not reported for the reason listed.								

Acoustic Output Reporting Table — Track 3, FDA 510(k) and IEC 60601-2-37

(Per transducer/mode that exceeds MI or TI value of 1.0)

Transducer model: **13L4**Operating mode: **M-mode**

Index Label			MI	TIS		TIB		TIC
				At surface	Below surface	At surface	Below surface	
Maximum Index Value			1.43	1.15		1.15		1.06
Index Component Value				1.15	0.62	1.15	0.75	
Acoustic Parameters	IEC	Units						
	$p_{ra,\alpha}$ at z_{MI}	(MPa)	3.69					
	P	(mW)		#		37.42		28.85
	P_{1x1}	(mW)		37.42		37.42		
	z_s	(cm)			1.40			
	z_b	(cm)					1.40	
	z_{MI}	(cm)	0.50					
	$z_{pii,\alpha}$	(cm)	0.70					
	f_{awf}	(MHz)	6.76	6.44		6.44		6.37
Other Information	p_{rr}	(Hz)	720.05					
	s_{rr}	(Hz)	#					
	η_{pps}	-	#					
	$I_{pa,\alpha}$ at $z_{pii,\alpha}$	(W/cm ²)	297.22					
	$I_{spta,\alpha}$ at $z_{pii,\alpha}$ or $z_{sii,\alpha}$	(mW/cm ²)	88.28					
	I_{spta} at z_{pii} or z_{sii}	(mW/cm ²)	122.34					
	p_r at z_{pii}	(MPa)	2.77					
	Operating Control Conditions	TX Level	(dB)	0.00	0.00		0.00	
Focus		(mm)	5.00	70.00		70.00		40.00
a This Index is not relevant to this operating mode.								
b This transducer is not intended for transcranial or neonatal cephalic uses.								
c This formulation for TIS is less than that for an alternate formulation in this mode.								
# No data is provided for this operation condition since the maximum index value is not reported for the reason listed.								

Acoustic Output Reporting Table — Track 3, FDA 510(k) and IEC 60601-2-37

(Per transducer/mode that exceeds MI or TI value of 1.0)

Transducer model: **13L4**Operating mode: **Pulsed Wave Doppler**

Index Label			MI	TIS		TIB		TIC
				At surface	Below surface	At surface	Below surface	
Maximum Index Value			1.43	1.67		1.73		2.17
Index Component Value				1.67	0.83	0.63	1.73	
Acoustic Parameters	IEC	Units						
	$p_{ra,\alpha}$ at z_{MI}	(MPa)	3.31					
	P	(mW)		#		24.53		24.53
	P_{1x1}	(mW)		58.92		24.53		
	z_s	(cm)			1.70			
	z_b	(cm)					0.50	
	z_{MI}	(cm)	0.50					
	$z_{pii,\alpha}$	(cm)	0.60					
	f_{awf}	(MHz)	5.35	5.94		5.39		5.39
Other Information	p_{rr}	(Hz)	151.99					
	s_{rr}	(Hz)	#					
	η_{pps}	-	#					
	$I_{pa,\alpha}$ at $z_{pii,\alpha}$	(W/cm ²)	456.45					
	$I_{spta,\alpha}$ at $z_{pii,\alpha}$ or $z_{sii,\alpha}$	(mW/cm ²)	59.11					
	I_{spta} at z_{pii} or z_{sii}	(mW/cm ²)	73.76					
	p_r at z_{pii}	(MPa)	3.55					
Operating Control Conditions	TX Level	(dB)	0.00	0.00		0.00		0.00
	Focus	(mm)	10.00	70.00		5.00		5.00
a	This Index is not relevant to this operating mode.							
b	This transducer is not intended for transcranial or neonatal cephalic uses.							
c	This formulation for TIS is less than that for an alternate formulation in this mode.							
#	No data is provided for this operation condition since the maximum index value is not reported for the reason listed.							

Acoustic Output Reporting Table — Track 3, FDA 510(k) and IEC 60601-2-37

(Per transducer/mode that exceeds MI or TI value of 1.0)

Transducer model: **13L4**Operating mode: **Color Flow Doppler**

Index Label			MI	TIS		TIB		TIC
				At surface	Below surface	At surface	Below surface	
Maximum Index Value			1.45	1.01		1.01		1.48
Index Component Value				1.01	1.01	1.01	1.01	
Acoustic Parameters	IEC	Units						
	$p_{ra,\alpha}$ at z_{MI}	(MPa)	3.35					
	P	(mW)		#		44.24		87.76
	$P_{1\times1}$	(mW)		36.21		36.21		
	z_s	(cm)			#			
	z_b	(cm)					#	
	z_{MI}	(cm)	1.30					
	$z_{pii,\alpha}$	(cm)	1.30					
	f_{awf}	(MHz)	5.27	5.88		5.88		5.94
Other Information	p_{rr}	(Hz)	#					
	s_{rr}	(Hz)	187.33					
	n_{pps}	-	12.00					
	$I_{pa,\alpha}$ at $z_{pii,\alpha}$	(W/cm ²)	520.98					
	$I_{spta,\alpha}$ at $z_{pii,\alpha}$ or $z_{sii,\alpha}$	(mW/cm ²)	21.21					
	I_{spta} at z_{pii} or z_{sii}	(mW/cm ²)	29.41					
	p_r at z_{pii}	(MPa)	4.17					
Operating Control Conditions	TX Level	(dB)	0.00	0.00		0.00		0.00
	Focus	(mm)	21.00	24.70		24.70		69.60
a	This Index is not relevant to this operating mode.							
b	This transducer is not intended for transcranial or neonatal cephalic uses.							
c	This formulation for TIS is less than that for an alternate formulation in this mode.							
#	No data is provided for this operation condition since the maximum index value is not reported for the reason listed.							

Acoustic Output Reporting Table — Track 3, FDA 510(k) and IEC 60601-2-37

(Per transducer/mode that exceeds MI or TI value of 1.0)

Transducer model: **14L4**Operating mode: **2D-mode**

Index Label			MI	TIS		TIB		TIC
				At surface	Below surface	At surface	Below surface	
Maximum Index Value			1.40	0.90		0.90		1.11
Index Component Value				0.90	0.90	0.90	0.90	
Acoustic Parameters	IEC	Units						
	$p_{ra,\alpha}$ at z_{MI}	(MPa)	3.22					
	P	(mW)		#		43.74		71.75
	P_{1x1}	(mW)		27.21		27.21		
	z_s	(cm)			#			
	z_b	(cm)					#	
	z_{MI}	(cm)	1.00					
	$z_{pii,\alpha}$	(cm)	1.00					
	f_{awf}	(MHz)	7.09	6.97		6.97		6.97
Other Information	p_{rr}	(Hz)	#					
	s_{rr}	(Hz)	52.13					
	n_{pps}	-	1.00					
	$I_{pa,\alpha}$ at $z_{pii,\alpha}$	(W/cm ²)	399.55					
	$I_{spta,\alpha}$ at $z_{pii,\alpha}$ or $z_{sii,\alpha}$	(mW/cm ²)	8.56					
	I_{spta} at z_{pii} or z_{sii}	(mW/cm ²)	12.05					
	p_r at z_{pii}	(MPa)	3.94					
Operating Control Conditions	TX Level	(dB)	0.00	0.00		0.00		0.00
	Focus	(mm)	15.00	60.00		60.00		60.00
a	This Index is not relevant to this operating mode.							
b	This transducer is not intended for transcranial or neonatal cephalic uses.							
c	This formulation for TIS is less than that for an alternate formulation in this mode.							
#	No data is provided for this operation condition since the maximum index value is not reported for the reason listed.							

Acoustic Output Reporting Table — Track 3, FDA 510(k) and IEC 60601-2-37

(Per transducer/mode that exceeds MI or TI value of 1.0)

Transducer model: **14L4**Operating mode: **M-mode**

Index Label			MI	TIS		TIB		TIC
				At surface	Below surface	At surface	Below surface	
Maximum Index Value			1.40	0.81		0.82		0.92
Index Component Value				0.81	0.49	0.65	0.82	
Acoustic Parameters	IEC	Units						
	$p_{ra,\alpha}$ at z_{MI}	(MPa)	3.19					
	P	(mW)		#		28.81		37.18
	$P_{1\times1}$	(mW)		37.18		#		
	z_s	(cm)			1.60			
	z_b	(cm)					1.40	
	z_{MI}	(cm)	1.20					
	$z_{pii,\alpha}$	(cm)	1.20					
	f_{awf}	(MHz)	6.94	4.55		4.72		4.55
Other Information	p_{rr}	(Hz)	720.00					
	s_{rr}	(Hz)	#					
	n_{pps}	-	#					
	$I_{pa,\alpha}$ at $z_{pii,\alpha}$	(W/cm ²)	400.04					
	$I_{spta,\alpha}$ at $z_{pii,\alpha}$ or $z_{sii,\alpha}$	(mW/cm ²)	62.80					
	I_{spta} at z_{pii} or z_{sii}	(mW/cm ²)	102.57					
	p_r at z_{pii}	(MPa)	3.88					
Operating Control Conditions	TX Level	(dB)	0.00	0.00		0.00		0.00
	Focus	(mm)	20.00	80.00		60.00		80.00
a	This Index is not relevant to this operating mode.							
b	This transducer is not intended for transcranial or neonatal cephalic uses.							
c	This formulation for TIS is less than that for an alternate formulation in this mode.							
#	No data is provided for this operation condition since the maximum index value is not reported for the reason listed.							

Acoustic Output Reporting Table — Track 3, FDA 510(k) and IEC 60601-2-37

(Per transducer/mode that exceeds MI or TI value of 1.0)

Transducer model: **14L4**Operating mode: **Pulsed Wave Doppler**

Index Label			MI	TIS		TIB		TIC
				At surface	Below surface	At surface	Below surface	
Maximum Index Value			1.39	1.23		1.55		1.36
Index Component Value				1.23	0.63	0.67	1.55	
Acoustic Parameters	IEC	Units						
	$p_{ra,\alpha}$ at z_{MI}	(MPa)	3.19					
	P	(mW)		#		26.52		13.48
	P_{1x1}	(mW)		42.67		26.52		
	z_s	(cm)			1.60			
	z_b	(cm)					1.10	
	z_{MI}	(cm)	1.20					
	$z_{pii,\alpha}$	(cm)	1.30					
	f_{awf}	(MHz)	5.36	6.05		5.32		5.39
Other Information	p_{rr}	(Hz)	243.00					
	s_{rr}	(Hz)	#					
	n_{pps}	-	#					
	$I_{pa,\alpha}$ at $z_{pii,\alpha}$	(W/cm ²)	405.36					
	$I_{spta,\alpha}$ at $z_{pii,\alpha}$ or $z_{sii,\alpha}$	(mW/cm ²)	74.05					
	I_{spta} at z_{pii} or z_{sii}	(mW/cm ²)	112.52					
	p_r at z_{pii}	(MPa)	3.68					
Operating Control Conditions	TX Level	(dB)	0.00	0.00		0.00		0.00
	Focus	(mm)	60.00	80.00		30.00		5.00
a This Index is not relevant to this operating mode.								
b This transducer is not intended for transcranial or neonatal cephalic uses.								
c This formulation for TIS is less than that for an alternate formulation in this mode.								
# No data is provided for this operation condition since the maximum index value is not reported for the reason listed.								

Acoustic Output Reporting Table — Track 3, FDA 510(k) and IEC 60601-2-37

(Per transducer/mode that exceeds MI or TI value of 1.0)

Transducer model: **14L4**Operating mode: **Color Flow Doppler**

Index Label			MI	TIS		TIB		TIC
				At surface	Below surface	At surface	Below surface	
Maximum Index Value			1.35	0.78		0.78		1.04
Index Component Value				0.77	0.77	0.77	0.77	
Acoustic Parameters	IEC	Units						
	$p_{ra,\alpha}$ at z_{MI}	(MPa)	2.86					
	P	(mW)		#		28.82		46.41
	$P_{1\times1}$	(mW)		26.44		26.44		
	z_s	(cm)			#			
	z_b	(cm)					#	
	z_{MI}	(cm)	0.90					
	$z_{pii,\alpha}$	(cm)	0.90					
	f_{awf}	(MHz)	4.51	6.10		6.10		5.37
Other Information	prr	(Hz)	#					
	srr	(Hz)	186.58					
	n_{pps}	-	12.00					
	$I_{pa,\alpha}$ at $z_{pii,\alpha}$	(W/cm ²)	406.85					
	$I_{spta,\alpha}$ at $z_{pii,\alpha}$ or $z_{sii,\alpha}$	(mW/cm ²)	4.72					
	I_{spta} at z_{pii} or z_{sii}	(mW/cm ²)	6.05					
	p_r at z_{pii}	(MPa)	3.22					
	Operating Control Conditions	TX Level	(dB)	0.00	0.00		0.00	
Focus		(mm)	15.00	10.00		10.00		10.00
a	This Index is not relevant to this operating mode.							
b	This transducer is not intended for transcranial or neonatal cephalic uses.							
c	This formulation for TIS is less than that for an alternate formulation in this mode.							
#	No data is provided for this operation condition since the maximum index value is not reported for the reason listed.							

Acoustic Output Reporting Table — Track 3, FDA 510(k) and IEC 60601-2-37

(Per transducer/mode that exceeds MI or TI value of 1.0)

Transducer model: **14L4**Operating mode: **Contrast**

Index Label			MI	TIS		TIB		TIC
				At surface	Below surface	At surface	Below surface	
Maximum Index Value			1.40	0.80		0.80		0.97
Index Component Value				0.52	0.52	0.52	0.52	
Acoustic Parameters	IEC	Units						
	$p_{ra,\alpha}$ at z_{MI}	(MPa)	3.19					
	P	(mW)		#		23.65		41.96
	P_{1x1}	(mW)		16.74		16.74		
	z_s	(cm)			#			
	z_b	(cm)					#	
	z_{MI}	(cm)	1.20					
	$z_{pii,\alpha}$	(cm)	1.20					
	f_{awf}	(MHz)	6.94	6.51		6.51		6.51
Other Information	p_{rr}	(Hz)	#					
	s_{rr}	(Hz)	28.53					
	n_{pps}	-	1.00					
	$I_{pa,\alpha}$ at $z_{pii,\alpha}$	(W/cm ²)	400.04					
	$I_{spta,\alpha}$ at $z_{pii,\alpha}$ or $z_{sii,\alpha}$	(mW/cm ²)	7.53					
	I_{spta} at z_{pii} or z_{sii}	(mW/cm ²)	11.04					
	p_r at z_{pii}	(MPa)	3.88					
Operating Control Conditions	TX Level	(dB)	0.00	0.00		0.00		0.00
	Focus	(mm)	20.00	20.00		20.00		60.00
a This Index is not relevant to this operating mode.								
b This transducer is not intended for transcranial or neonatal cephalic uses.								
c This formulation for TIS is less than that for an alternate formulation in this mode.								
# No data is provided for this operation condition since the maximum index value is not reported for the reason listed.								

Acoustic Output Reporting Table — Track 3, FDA 510(k) and IEC 60601-2-37

(Per transducer/mode that exceeds MI or TI value of 1.0)

Transducer model: **16L4**Operating mode: **2D-mode**

Index Label			MI	TIS		TIB		TIC
				At surface	Below surface	At surface	Below surface	
Maximum Index Value			1.20	1.03		1.03		1.20
Index Component Value				1.03	1.03	1.03	1.03	
Acoustic Parameters	IEC	Units						
	$p_{ra,\alpha}$ at z_{MI}	(MPa)	3.06					
	P	(mW)		#		34.14		51.44
	$P_{1\times1}$	(mW)		26.49		26.49		
	z_s	(cm)			#			
	z_b	(cm)					#	
	z_{MI}	(cm)	1.00					
	$z_{pii,\alpha}$	(cm)	1.00					
	f_{awf}	(MHz)	6.84	8.19		8.19		8.09
Other Information	p_{rr}	(Hz)	#					
	s_{rr}	(Hz)	116.60					
	n_{pps}	-	2.00					
	$I_{pa,\alpha}$ at $z_{pii,\alpha}$	(W/cm ²)	428.79					
	$I_{spta,\alpha}$ at $z_{pii,\alpha}$ or $z_{sii,\alpha}$	(mW/cm ²)	19.63					
	I_{spta} at z_{pii} or z_{sii}	(mW/cm ²)	27.31					
	p_r at z_{pii}	(MPa)	3.63					
Operating Control Conditions	TX Level	(dB)	0.00	0.00		0.00		0.00
	Focus	(mm)	12.50	30.00		30.00		40.00
a	This Index is not relevant to this operating mode.							
b	This transducer is not intended for transcranial or neonatal cephalic uses.							
c	This formulation for TIS is less than that for an alternate formulation in this mode.							
#	No data is provided for this operation condition since the maximum index value is not reported for the reason listed.							

Acoustic Output Reporting Table — Track 3, FDA 510(k) and IEC 60601-2-37

(Per transducer/mode that exceeds MI or TI value of 1.0)

Transducer model: **16L4**Operating mode: **M-mode**

Index Label			MI	TIS		TIB		TIC
				At surface	Below surface	At surface	Below surface	
Maximum Index Value			1.05	0.49		0.49		0.52
Index Component Value				0.49	0.28	0.49	0.36	
Acoustic Parameters	IEC	Units						
	$p_{ra,\alpha}$ at z_{MI}	(MPa)	2.95					
	P	(mW)		#		12.74		12.74
	P_{1x1}	(mW)		12.74		12.74		
	z_s	(cm)			1.00			
	z_b	(cm)					1.00	
	z_{MI}	(cm)	0.70					
	$z_{pii,\alpha}$	(cm)	0.80					
	f_{awf}	(MHz)	8.19	8.09		8.09		8.09
Other Information	p_{rr}	(Hz)	720.00					
	s_{rr}	(Hz)	#					
	n_{pps}	-	#					
	$I_{pa,\alpha}$ at $z_{pii,\alpha}$	(W/cm ²)	259.44					
	$I_{spta,\alpha}$ at $z_{pii,\alpha}$ or $z_{sii,\alpha}$	(mW/cm ²)	50.14					
	I_{spta} at z_{pii} or z_{sii}	(mW/cm ²)	78.81					
	p_r at z_{pii}	(MPa)	3.37					
Operating Control Conditions	TX Level	(dB)	0.00	0.00		0.00		0.00
	Focus	(mm)	30.00	40.00		40.00		40.00
a This Index is not relevant to this operating mode.								
b This transducer is not intended for transcranial or neonatal cephalic uses.								
c This formulation for TIS is less than that for an alternate formulation in this mode.								
# No data is provided for this operation condition since the maximum index value is not reported for the reason listed.								

Acoustic Output Reporting Table — Track 3, FDA 510(k) and IEC 60601-2-37

(Per transducer/mode that exceeds MI or TI value of 1.0)

Transducer model: **16L4**Operating mode: **Pulsed Wave Doppler**

Index Label			MI	TIS		TIB		TIC
				At surface	Below surface	At surface	Below surface	
Maximum Index Value			1.48	0.84		1.92		1.26
Index Component Value				0.84	0.55	0.26	1.92	
Acoustic Parameters	IEC	Units						
	$p_{ra,\alpha}$ at z_{MI}	(MPa)	3.67					
	P	(mW)		#		21.51		22.11
	$P_{1\times1}$	(mW)		28.72		21.51		
	z_s	(cm)			1.00			
	z_b	(cm)					0.40	
	z_{MI}	(cm)	0.80					
	$z_{pii,\alpha}$	(cm)	0.80					
	f_{awf}	(MHz)	6.16	6.15		5.52		6.13
Other Information	prr	(Hz)	976.56					
	srr	(Hz)	#					
	n_{pps}	-	#					
	$I_{pa,\alpha}$ at $z_{pii,\alpha}$	(W/cm ²)	703.10					
	$I_{spta,\alpha}$ at $z_{pii,\alpha}$ or $z_{sii,\alpha}$	(mW/cm ²)	426.54					
	I_{spta} at z_{pii} or z_{sii}	(mW/cm ²)	581.72					
	p_r at z_{pii}	(MPa)	4.10					
Operating Control Conditions	TX Level	(dB)	0.00	0.00		0.00		0.00
	Focus	(mm)	12.50	40.00		5.00		20.00
a	This Index is not relevant to this operating mode.							
b	This transducer is not intended for transcranial or neonatal cephalic uses.							
c	This formulation for TIS is less than that for an alternate formulation in this mode.							
#	No data is provided for this operation condition since the maximum index value is not reported for the reason listed.							

Acoustic Output Reporting Table — Track 3, FDA 510(k) and IEC 60601-2-37

(Per transducer/mode that exceeds MI or TI value of 1.0)

Transducer model: **16L4**Operating mode: **Color Flow Doppler**

Index Label			MI	TIS		TIB		TIC
				At surface	Below surface	At surface	Below surface	
Maximum Index Value			1.50	0.70		0.70		1.07
Index Component Value				0.70	0.70	0.70	0.70	
Acoustic Parameters	IEC	Units						
	$p_{ra,\alpha}$ at z_{MI}	(MPa)	3.70					
	P	(mW)		#		27.11		39.52
	P_{1x1}	(mW)		23.91		23.91		
	z_s	(cm)			#			
	z_b	(cm)					#	
	z_{MI}	(cm)	1.00					
	$z_{pii,\alpha}$	(cm)	1.00					
	f_{awf}	(MHz)	6.15	6.17		6.17		6.18
Other Information	p_{rr}	(Hz)	#					
	s_{rr}	(Hz)	85.60					
	n_{pps}	-	12.00					
	$I_{pa,\alpha}$ at $z_{pii,\alpha}$	(W/cm ²)	760.83					
	$I_{spta,\alpha}$ at $z_{pii,\alpha}$ or $z_{sii,\alpha}$	(mW/cm ²)	6.09					
	I_{spta} at z_{pii} or z_{sii}	(mW/cm ²)	8.20					
	p_r at z_{pii}	(MPa)	4.58					
Operating Control Conditions	TX Level	(dB)	0.00	0.00		0.00		0.00
	Focus	(mm)	15.00	7.50		7.50		5.00
a This Index is not relevant to this operating mode.								
b This transducer is not intended for transcranial or neonatal cephalic uses.								
c This formulation for TIS is less than that for an alternate formulation in this mode.								
# No data is provided for this operation condition since the maximum index value is not reported for the reason listed.								

Acoustic Output Reporting Table — Track 3, FDA 510(k) and IEC 60601-2-37

(Per transducer/mode that exceeds MI or TI value of 1.0)

Transducer model: **18H5**Operating mode: **2D-mode**

Index Label			MI	TIS		TIB		TIC
				At surface	Below surface	At surface	Below surface	
Maximum Index Value			1.49	1.56		1.56		3.14
Index Component Value				1.56	1.56	1.56	1.56	
Acoustic Parameters	IEC	Units						
	$p_{ra,\alpha}$ at z_{MI}	(MPa)	4.02					
	P	(mW)		#		118.63		118.63
	$P_{1\times1}$	(mW)		44.34		44.34		
	z_s	(cm)			#			
	z_b	(cm)					#	
	z_{MI}	(cm)	0.70					
	$z_{pii,\alpha}$	(cm)	0.70					
	f_{awf}	(MHz)	7.58	7.38		7.38		7.38
Other Information	p_{rr}	(Hz)	#					
	s_{rr}	(Hz)	38.02					
	n_{pps}	-	1.00					
	$I_{pa,\alpha}$ at $z_{pii,\alpha}$	(W/cm ²)	599.81					
	$I_{spta,\alpha}$ at $z_{pii,\alpha}$ or $z_{sii,\alpha}$	(mW/cm ²)	21.58					
	I_{spta} at z_{pii} or z_{sii}	(mW/cm ²)	31.13					
	p_r at z_{pii}	(MPa)	4.57					
Operating Control Conditions	TX Level	(dB)	0.00	0.00		0.00		0.00
	Focus	(mm)	15.00	60.00		60.00		60.00
a This Index is not relevant to this operating mode.								
b This transducer is not intended for transcranial or neonatal cephalic uses.								
c This formulation for TIS is less than that for an alternate formulation in this mode.								
# No data is provided for this operation condition since the maximum index value is not reported for the reason listed.								

Acoustic Output Reporting Table — Track 3, FDA 510(k) and IEC 60601-2-37

(Per transducer/mode that exceeds MI or TI value of 1.0)

Transducer model: **18H5**Operating mode: **M-mode**

Index Label			MI	TIS		TIB		TIC
				At surface	Below surface	At surface	Below surface	
Maximum Index Value			1.49	0.86		0.86		0.94
Index Component Value				0.86	0.52	0.86	0.48	
Acoustic Parameters	IEC	Units						
	$p_{ra,\alpha}$ at z_{MI}	(MPa)	4.02					
	P	(mW)		#		24.57		24.57
	P_{1x1}	(mW)		24.57		24.57		
	z_s	(cm)			1.00			
	z_b	(cm)					1.00	
	z_{MI}	(cm)	0.70					
	$z_{pii,\alpha}$	(cm)	0.70					
	f_{awf}	(MHz)	7.58	7.38		7.38		7.38
Other Information	p_{rr}	(Hz)	720.00					
	s_{rr}	(Hz)	#					
	n_{pps}	-	#					
	$I_{pa,\alpha}$ at $z_{pii,\alpha}$	(W/cm ²)	599.81					
	$I_{spta,\alpha}$ at $z_{pii,\alpha}$ or $z_{sii,\alpha}$	(mW/cm ²)	81.56					
	I_{spta} at z_{pii} or z_{sii}	(mW/cm ²)	115.64					
	p_r at z_{pii}	(MPa)	4.57					
Operating Control Conditions	TX Level	(dB)	0.00	0.00		0.00		0.00
	Focus	(mm)	15.00	60.00		60.00		60.00
a This Index is not relevant to this operating mode.								
b This transducer is not intended for transcranial or neonatal cephalic uses.								
c This formulation for TIS is less than that for an alternate formulation in this mode.								
# No data is provided for this operation condition since the maximum index value is not reported for the reason listed.								

Acoustic Output Reporting Table — Track 3, FDA 510(k) and IEC 60601-2-37

(Per transducer/mode that exceeds MI or TI value of 1.0)

Transducer model: **18H5**Operating mode: **Pulsed Wave Doppler**

Index Label			MI	TIS		TIB		TIC
				At surface	Below surface	At surface	Below surface	
Maximum Index Value			1.50	2.15		2.15		2.15
Index Component Value				2.15	1.20	2.15	1.33	
Acoustic Parameters	IEC	Units						
	$p_{ra,\alpha}$ at z_{MI}	(MPa)	3.97					
	P	(mW)		#		64.40		45.69
	$P_{1\times1}$	(mW)		64.40		64.40		
	z_s	(cm)			1.20			
	z_b	(cm)					1.20	
	z_{MI}	(cm)	0.80					
	$z_{pii,\alpha}$	(cm)	0.80					
	f_{awf}	(MHz)	7.00	7.00		7.00		7.04
Other Information	prr	(Hz)	1220.70					
	srr	(Hz)	#					
	n_{pps}	-	#					
	$I_{pa,\alpha}$ at $z_{pii,\alpha}$	(W/cm ²)	662.70					
	$I_{spta,\alpha}$ at $z_{pii,\alpha}$ or $z_{sii,\alpha}$	(mW/cm ²)	480.28					
	I_{spta} at z_{pii} or z_{sii}	(mW/cm ²)	706.86					
	p_r at z_{pii}	(MPa)	4.82					
Operating Control Conditions	TX Level	(dB)	0.00	0.00		0.00		0.00
	Focus	(mm)	15.00	60.00		60.00		30.00
a	This Index is not relevant to this operating mode.							
b	This transducer is not intended for transcranial or neonatal cephalic uses.							
c	This formulation for TIS is less than that for an alternate formulation in this mode.							
#	No data is provided for this operation condition since the maximum index value is not reported for the reason listed.							

Acoustic Output Reporting Table — Track 3, FDA 510(k) and IEC 60601-2-37

(Per transducer/mode that exceeds MI or TI value of 1.0)

Transducer model: **18H5**Operating mode: **Color Flow Doppler**

Index Label			MI	TIS		TIB		TIC
				At surface	Below surface	At surface	Below surface	
Maximum Index Value			1.50	1.40		1.40		1.24
Index Component Value				1.40	1.40	1.40	1.40	
Acoustic Parameters	IEC	Units						
	$p_{ra,\alpha}$ at z_{MI}	(MPa)	4.00					
	P	(mW)		#		43.94		46.84
	P_{1x1}	(mW)		41.44		41.44		
	z_s	(cm)			#			
	z_b	(cm)					#	
	z_{MI}	(cm)	0.60					
	$z_{pii,\alpha}$	(cm)	0.60					
	f_{awf}	(MHz)	7.06	7.10		7.10		7.12
Other Information	p_{rr}	(Hz)	#					
	s_{rr}	(Hz)	550.06					
	n_{pps}	-	16.00					
	$I_{pa,\alpha}$ at $z_{pii,\alpha}$	(W/cm ²)	799.35					
	$I_{spta,\alpha}$ at $z_{pii,\alpha}$ or $z_{sii,\alpha}$	(mW/cm ²)	148.08					
	I_{spta} at z_{pii} or z_{sii}	(mW/cm ²)	198.37					
	p_r at z_{pii}	(MPa)	4.63					
Operating Control Conditions	TX Level	(dB)	0.00	0.00		0.00		0.00
	Focus	(mm)	10.00	7.50		7.50		5.00
a This Index is not relevant to this operating mode.								
b This transducer is not intended for transcranial or neonatal cephalic uses.								
c This formulation for TIS is less than that for an alternate formulation in this mode.								
# No data is provided for this operation condition since the maximum index value is not reported for the reason listed.								

Acoustic Output Reporting Table — Track 3, FDA 510(k) and IEC 60601-2-37

(Per transducer/mode that exceeds MI or TI value of 1.0)

Transducer model: **11M3**Operating mode: **2D-mode**

Index Label			MI	TIS		TIB		TIC
				At surface	Below surface	At surface	Below surface	
Maximum Index Value			1.50	1.93		1.93		3.11
Index Component Value				1.93	1.93	1.93	1.93	
Acoustic Parameters	IEC	Units						
	$p_{ra,\alpha}$ at z_{MI}	(MPa)	3.18					
	P	(mW)		#		102.72		163.48
	$P_{1\times1}$	(mW)		72.95		72.95		
	z_s	(cm)			#			
	z_b	(cm)					#	
	z_{MI}	(cm)	1.40					
	$z_{pII,\alpha}$	(cm)	2.30					
	f_{awf}	(MHz)	4.52	5.55		5.55		5.55
Other Information	prr	(Hz)	#					
	srr	(Hz)	79.97					
	n_{pps}	-	1.00					
	$I_{pa,\alpha}$ at $z_{pII,\alpha}$	(W/cm ²)	155.98					
	$I_{spta,\alpha}$ at $z_{pII,\alpha}$ or $z_{sII,\alpha}$	(mW/cm ²)	57.09					
	I_{spta} at z_{pII} or z_{sII}	(mW/cm ²)	117.05					
	p_r at z_{pII}	(MPa)	2.75					
Operating Control Conditions	TX Level	(dB)	0.00	0.00		0.00		0.00
	Focus	(mm)	100.00	40.00		40.00		60.00
a	This Index is not relevant to this operating mode.							
b	This transducer is not intended for transcranial or neonatal cephalic uses.							
c	This formulation for TIS is less than that for an alternate formulation in this mode.							
#	No data is provided for this operation condition since the maximum index value is not reported for the reason listed.							

Acoustic Output Reporting Table — Track 3, FDA 510(k) and IEC 60601-2-37

(Per transducer/mode that exceeds MI or TI value of 1.0)

Transducer model: **11M3**Operating mode: **M-mode**

Index Label			MI	TIS		TIB		TIC
				At surface	Below surface	At surface	Below surface	
Maximum Index Value			1.50	2.58		2.58		2.38
Index Component Value				2.58	1.07	2.58	0.85	
Acoustic Parameters	IEC	Units						
	$p_{ra,\alpha}$ at z_{MI}	(MPa)	3.19					
	P	(mW)		#		98.12		91.39
	P_{1x1}	(mW)		98.12		#		
	z_s	(cm)			2.30			
	z_b	(cm)					2.30	
	z_{MI}	(cm)	1.40					
	$z_{pii,\alpha}$	(cm)	2.30					
	f_{awf}	(MHz)	4.53	5.52		5.52		4.38
Other Information	p_{rr}	(Hz)	720.00					
	s_{rr}	(Hz)	#					
	n_{pps}	-	#					
	$I_{pa,\alpha}$ at $z_{pii,\alpha}$	(W/cm ²)	158.57					
	$I_{spta,\alpha}$ at $z_{pii,\alpha}$ or $z_{sii,\alpha}$	(mW/cm ²)	79.58					
	I_{spta} at z_{pii} or z_{sii}	(mW/cm ²)	163.26					
	p_r at z_{pii}	(MPa)	2.78					
Operating Control Conditions	TX Level	(dB)	0.00	0.00		0.00		0.00
	Focus	(mm)	80.00	80.00		80.00		60.00
a This Index is not relevant to this operating mode.								
b This transducer is not intended for transcranial or neonatal cephalic uses.								
c This formulation for TIS is less than that for an alternate formulation in this mode.								
# No data is provided for this operation condition since the maximum index value is not reported for the reason listed.								

Acoustic Output Reporting Table — Track 3, FDA 510(k) and IEC 60601-2-37

(Per transducer/mode that exceeds MI or TI value of 1.0)

Transducer model: **11M3**Operating mode: **Pulsed Wave Doppler**

Index Label			MI	TIS		TIB		TIC
				At surface	Below surface	At surface	Below surface	
Maximum Index Value			1.46	1.22		1.84		1.81
Index Component Value				1.22	0.63	0.73	1.84	
Acoustic Parameters	IEC	Units						
	$p_{ra,\alpha}$ at z_{MI}	(MPa)	3.09					
	P	(mW)		#		34.87		16.78
	$P_{1\times1}$	(mW)		58.43		#		
	z_s	(cm)			2.50			
	z_b	(cm)					1.10	
	z_{MI}	(cm)	1.20					
	$z_{pii,\alpha}$	(cm)	1.20					
	f_{awf}	(MHz)	4.42	4.38		4.39		4.46
Other Information	pr_r	(Hz)	976.56					
	sr_r	(Hz)	#					
	n_{pps}	-	#					
	$I_{pa,\alpha}$ at $z_{pii,\alpha}$	(W/cm ²)	502.59					
	$I_{spta,\alpha}$ at $z_{pii,\alpha}$ or $z_{sii,\alpha}$	(mW/cm ²)	402.82					
	I_{spta} at z_{pii} or z_{sii}	(mW/cm ²)	580.99					
	p_r at z_{pii}	(MPa)	3.71					
Operating Control Conditions	TX Level	(dB)	0.00	0.00		0.00		0.00
	Focus	(mm)	20.00	80.00		40.00		5.00
a This Index is not relevant to this operating mode.								
b This transducer is not intended for transcranial or neonatal cephalic uses.								
c This formulation for TIS is less than that for an alternate formulation in this mode.								
# No data is provided for this operation condition since the maximum index value is not reported for the reason listed.								

Acoustic Output Reporting Table — Track 3, FDA 510(k) and IEC 60601-2-37

(Per transducer/mode that exceeds MI or TI value of 1.0)

Transducer model: **11M3**Operating mode: **Color Flow Doppler**

Index Label			MI	TIS		TIB		TIC
				At surface	Below surface	At surface	Below surface	
Maximum Index Value			1.50	1.20		1.20		2.27
Index Component Value				0.41	0.41	0.41	0.41	
Acoustic Parameters	IEC	Units						
	$p_{ra,\alpha}$ at z_{MI}	(MPa)	3.16					
	P	(mW)		#		20.94		37.64
	$P_{1\times1}$	(mW)		19.34		19.34		
	z_s	(cm)			#			
	z_b	(cm)					#	
	z_{MI}	(cm)	1.10					
	$z_{pii,\alpha}$	(cm)	1.10					
	f_{awf}	(MHz)	4.42	4.48		4.48		4.41
Other Information	p_{rr}	(Hz)	#					
	s_{rr}	(Hz)	325.35					
	n_{pps}	-	10.00					
	$I_{pa,\alpha}$ at $z_{pii,\alpha}$	(W/cm ²)	507.50					
	$I_{spta,\alpha}$ at $z_{pii,\alpha}$ or $z_{sii,\alpha}$	(mW/cm ²)	108.13					
	I_{spta} at z_{pii} or z_{sii}	(mW/cm ²)	151.19					
	p_r at z_{pii}	(MPa)	3.74					
Operating Control Conditions	TX Level	(dB)	0.00	0.00		0.00		0.00
	Focus	(mm)	20.00	40.00		40.00		120.00
a	This Index is not relevant to this operating mode.							
b	This transducer is not intended for transcranial or neonatal cephalic uses.							
c	This formulation for TIS is less than that for an alternate formulation in this mode.							
#	No data is provided for this operation condition since the maximum index value is not reported for the reason listed.							

Acoustic Output Reporting Table — Track 3, FDA 510(k) and IEC 60601-2-37

(Per transducer/mode that exceeds MI or TI value of 1.0)

Transducer model: **L10-5v**Operating mode: **2D-mode**

Index Label			MI	TIS		TIB		TIC
				At surface	Below surface	At surface	Below surface	
Maximum Index Value			1.36	1.07		1.07		1.12
Index Component Value				1.07	1.07	1.07	1.07	
Acoustic Parameters	IEC	Units						
	$p_{ra,\alpha}$ at z_{MI}	(MPa)	3.48					
	P	(mW)		#		45.61		70.19
	$P_{1\times1}$	(mW)		34.31		34.31		
	z_s	(cm)			#			
	z_b	(cm)					#	
	z_{MI}	(cm)	1.10					
	$z_{pii,\alpha}$	(cm)	1.10					
	f_{awf}	(MHz)	6.76	6.55		6.55		6.50
Other Information	prr	(Hz)	#					
	srr	(Hz)	135.87					
	n_{pps}	-	1.00					
	$I_{pa,\alpha}$ at $z_{pii,\alpha}$	(W/cm ²)	476.44					
	$I_{spta,\alpha}$ at $z_{pii,\alpha}$ or $z_{sii,\alpha}$	(mW/cm ²)	23.23					
	I_{spta} at z_{pii} or z_{sii}	(mW/cm ²)	33.74					
	p_r at z_{pii}	(MPa)	4.50					
Operating Control Conditions	TX Level	(dB)	0.00	0.00		0.00		0.00
	Focus	(mm)	15.00	50.00		50.00		70.00
a	This Index is not relevant to this operating mode.							
b	This transducer is not intended for transcranial or neonatal cephalic uses.							
c	This formulation for TIS is less than that for an alternate formulation in this mode.							
#	No data is provided for this operation condition since the maximum index value is not reported for the reason listed.							

Acoustic Output Reporting Table — Track 3, FDA 510(k) and IEC 60601-2-37

(Per transducer/mode that exceeds MI or TI value of 1.0)

Transducer model: **L10-5v**Operating mode: **M-mode**

Index Label			MI	TIS		TIB		TIC
				At surface	Below surface	At surface	Below surface	
Maximum Index Value			1.36	0.71		0.71		0.71
Index Component Value				0.71	0.40	0.71	0.71	
Acoustic Parameters	IEC	Units						
	$p_{ra,\alpha}$ at z_{MI}	(MPa)	3.48					
	P	(mW)		#		26.92		26.92
	P_{1x1}	(mW)		26.92		26.92		
	z_s	(cm)			1.50			
	z_b	(cm)					1.50	
	z_{MI}	(cm)	1.10					
	$z_{pii,\alpha}$	(cm)	1.10					
	f_{awf}	(MHz)	6.76	5.50		5.50		5.50
Other Information	p_{rr}	(Hz)	720.00					
	s_{rr}	(Hz)	#					
	n_{pps}	-	#					
	$I_{pa,\alpha}$ at $z_{pii,\alpha}$	(W/cm ²)	476.44					
	$I_{spta,\alpha}$ at $z_{pii,\alpha}$ or $z_{sii,\alpha}$	(mW/cm ²)	77.54					
	I_{spta} at z_{pii} or z_{sii}	(mW/cm ²)	129.55					
	p_r at z_{pii}	(MPa)	4.50					
Operating Control Conditions	TX Level	(dB)	0.00	0.00		0.00		0.00
	Focus	(mm)	15.00	70.00		70.00		70.00
a This Index is not relevant to this operating mode.								
b This transducer is not intended for transcranial or neonatal cephalic uses.								
c This formulation for TIS is less than that for an alternate formulation in this mode.								
# No data is provided for this operation condition since the maximum index value is not reported for the reason listed.								

Acoustic Output Reporting Table — Track 3, FDA 510(k) and IEC 60601-2-37

(Per transducer/mode that exceeds MI or TI value of 1.0)

Transducer model: **L10-5v**Operating mode: **Pulsed Wave Doppler**

Index Label			MI	TIS		TIB		TIC
				At surface	Below surface	At surface	Below surface	
Maximum Index Value			1.36	1.93		2.10		1.65
Index Component Value				1.93	0.99	1.45	2.10	
Acoustic Parameters	IEC	Units						
	$p_{ra,\alpha}$ at z_{MI}	(MPa)	2.74					
	P	(mW)		#		51.17		60.33
	$P_{1\times1}$	(mW)		67.54		51.17		
	z_s	(cm)			1.60			
	z_b	(cm)					1.20	
	z_{MI}	(cm)	1.50					
	$z_{pii,\alpha}$	(cm)	1.50					
	f_{awf}	(MHz)	4.08	5.99		5.94		5.90
Other Information	pr_r	(Hz)	1562.50					
	sr_r	(Hz)	#					
	n_{pps}	-	#					
	$I_{pa,\alpha}$ at $z_{pii,\alpha}$	(W/cm ²)	250.07					
	$I_{spta,\alpha}$ at $z_{pii,\alpha}$ or $z_{sii,\alpha}$	(mW/cm ²)	294.80					
	I_{spta} at z_{pii} or z_{sii}	(mW/cm ²)	456.13					
	p_r at z_{pii}	(MPa)	3.30					
Operating Control Conditions	TX Level	(dB)	0.00	0.00		0.00		0.00
	Focus	(mm)	20.00	50.00		30.00		40.00
a	This Index is not relevant to this operating mode.							
b	This transducer is not intended for transcranial or neonatal cephalic uses.							
c	This formulation for TIS is less than that for an alternate formulation in this mode.							
#	No data is provided for this operation condition since the maximum index value is not reported for the reason listed.							

Acoustic Output Reporting Table — Track 3, FDA 510(k) and IEC 60601-2-37

(Per transducer/mode that exceeds MI or TI value of 1.0)

Transducer model: **L10-5v**Operating mode: **Color Flow Doppler**

Index Label			MI	TIS		TIB		TIC
				At surface	Below surface	At surface	Below surface	
Maximum Index Value			1.36	1.29		1.29		1.57
Index Component Value				1.21	1.21	1.21	1.21	
Acoustic Parameters	IEC	Units						
	$p_{ra,\alpha}$ at z_{MI}	(MPa)	3.18					
	P	(mW)		#		45.83		97.81
	P_{1x1}	(mW)		41.96		41.96		
	z_s	(cm)			#			
	z_b	(cm)					#	
	z_{MI}	(cm)	1.10					
	$z_{pii,\alpha}$	(cm)	1.10					
	f_{awf}	(MHz)	5.45	6.03		6.03		6.02
Other Information	p_{rr}	(Hz)	#					
	s_{rr}	(Hz)	185.74					
	n_{pps}	-	12.00					
	$I_{pa,\alpha}$ at $z_{pii,\alpha}$	(W/cm ²)	512.28					
	$I_{spta,\alpha}$ at $z_{pii,\alpha}$ or $z_{sii,\alpha}$	(mW/cm ²)	6.40					
	I_{spta} at z_{pii} or z_{sii}	(mW/cm ²)	9.32					
	p_r at z_{pii}	(MPa)	3.81					
Operating Control Conditions	TX Level	(dB)	0.00	0.00		0.00		0.00
	Focus	(mm)	15.00	7.50		7.50		70.00
a	This Index is not relevant to this operating mode.							
b	This transducer is not intended for transcranial or neonatal cephalic uses.							
c	This formulation for TIS is less than that for an alternate formulation in this mode.							
#	No data is provided for this operation condition since the maximum index value is not reported for the reason listed.							

Acoustic Output Reporting Table — Track 3, FDA 510(k) and IEC 60601-2-37

(Per transducer/mode that exceeds MI or TI value of 1.0)

Transducer model: **5P1**Operating mode: **2D-mode**

Index Label			MI	TIS		TIB		TIC
				At surface	Below surface	At surface	Below surface	
Maximum Index Value			1.50	0.91		0.91		1.65
Index Component Value				0.91	0.91	0.91	0.91	
Acoustic Parameters	IEC	Units						
	$p_{ra,\alpha}$ at z_{MI}	(MPa)	1.89					
	P	(mW)		#		75.52		117.84
	$P_{1\times1}$	(mW)		55.94		55.94		
	z_s	(cm)			#			
	z_b	(cm)					#	
	z_{MI}	(cm)	4.50					
	$z_{pii,\alpha}$	(cm)	4.50					
	f_{awf}	(MHz)	1.60	3.40		3.40		1.83
Other Information	p_{rr}	(Hz)	#					
	s_{rr}	(Hz)	299.76					
	n_{pps}	-	2.00					
	$I_{pa,\alpha}$ at $z_{pii,\alpha}$	(W/cm ²)	135.44					
	$I_{spta,\alpha}$ at $z_{pii,\alpha}$ or $z_{sii,\alpha}$	(mW/cm ²)	122.79					
	I_{spta} at z_{pii} or z_{sii}	(mW/cm ²)	174.98					
	p_r at z_{pii}	(MPa)	2.41					
Operating Control Conditions	TX Level	(dB)	0.00	0.00		0.00		0.00
	Focus	(mm)	80.00	40.00		40.00		80.00
a This Index is not relevant to this operating mode.								
b This transducer is not intended for transcranial or neonatal cephalic uses.								
c This formulation for TIS is less than that for an alternate formulation in this mode.								
# No data is provided for this operation condition since the maximum index value is not reported for the reason listed.								

Acoustic Output Reporting Table — Track 3, FDA 510(k) and IEC 60601-2-37

(Per transducer/mode that exceeds MI or TI value of 1.0)

Transducer model: **5P1**Operating mode: **M-mode**

Index Label			MI	TIS		TIB		TIC
				At surface	Below surface	At surface	Below surface	
Maximum Index Value			1.50	1.07		3.21		2.96
Index Component Value				0.74	1.07	0.63	3.21	
Acoustic Parameters	IEC	Units						
	$p_{ra,\alpha}$ at z_{MI}	(MPa)	2.01					
	P	(mW)		#		179.15		210.67
	P_{1x1}	(mW)		84.81		#		
	z_s	(cm)			4.30			
	z_b	(cm)					5.50	
	z_{MI}	(cm)	5.70					
	$z_{pii,\alpha}$	(cm)	5.70					
	f_{awf}	(MHz)	1.83	1.82		1.83		1.82
Other Information	p_{rr}	(Hz)	720.00					
	s_{rr}	(Hz)	#					
	n_{pps}	-	#					
	$I_{pa,\alpha}$ at $z_{pii,\alpha}$	(W/cm ²)	294.19					
	$I_{spta,\alpha}$ at $z_{pii,\alpha}$ or $z_{sii,\alpha}$	(mW/cm ²)	296.33					
	I_{spta} at z_{pii} or z_{sii}	(mW/cm ²)	587.12					
	p_r at z_{pii}	(MPa)	2.71					
Operating Control Conditions	TX Level	(dB)	0.00	0.00		0.00		0.00
	Focus	(mm)	80.00	100.00		80.00		100.00
a This Index is not relevant to this operating mode.								
b This transducer is not intended for transcranial or neonatal cephalic uses.								
c This formulation for TIS is less than that for an alternate formulation in this mode.								
# No data is provided for this operation condition since the maximum index value is not reported for the reason listed.								

Acoustic Output Reporting Table — Track 3, FDA 510(k) and IEC 60601-2-37

(Per transducer/mode that exceeds MI or TI value of 1.0)

Transducer model: **5P1**Operating mode: **Pulsed Wave Doppler**

Index Label			MI	TIS		TIB		TIC
				At surface	Below surface	At surface	Below surface	
Maximum Index Value			1.50	1.50		3.66		2.59
Index Component Value				1.50	1.24	0.39	3.66	
Acoustic Parameters	IEC	Units						
	$p_{ra,\alpha}$ at z_{MI}	(MPa)	2.35					
	P	(mW)		#		104.35		96.52
	$P_{1\times1}$	(mW)		90.93		#		
	z_s	(cm)			2.00			
	z_b	(cm)					3.90	
	z_{MI}	(cm)	4.50					
	$z_{pII,\alpha}$	(cm)	4.50					
	f_{awf}	(MHz)	2.47	3.46		1.57		3.49
Other Information	prr	(Hz)	610.35					
	srr	(Hz)	#					
	n_{pps}	-	#					
	$I_{pa,\alpha}$ at $z_{pII,\alpha}$	(W/cm ²)	448.17					
	$I_{spta,\alpha}$ at $z_{pII,\alpha}$ or $z_{sII,\alpha}$	(mW/cm ²)	447.63					
	I_{spta} at z_{pII} or z_{sII}	(mW/cm ²)	943.81					
	p_r at z_{pII}	(MPa)	3.37					
Operating Control Conditions	TX Level	(dB)	0.00	0.00		0.00		0.00
	Focus	(mm)	60.00	40.00		60.00		20.00
a	This Index is not relevant to this operating mode.							
b	This transducer is not intended for transcranial or neonatal cephalic uses.							
c	This formulation for TIS is less than that for an alternate formulation in this mode.							
#	No data is provided for this operation condition since the maximum index value is not reported for the reason listed.							

Acoustic Output Reporting Table — Track 3, FDA 510(k) and IEC 60601-2-37

(Per transducer/mode that exceeds MI or TI value of 1.0)

Transducer model: **5P1**Operating mode: **Color Flow Doppler**

Index Label			MI	TIS		TIB		TIC
				At surface	Below surface	At surface	Below surface	
Maximum Index Value			1.50	1.22		1.22		2.43
Index Component Value				0.82	0.82	0.82	0.82	
Acoustic Parameters	IEC	Units						
	$p_{ra,\alpha}$ at z_{MI}	(MPa)	2.35					
	P	(mW)		#		66.83		55.81
	P_{1x1}	(mW)		49.50		49.50		
	z_s	(cm)			#			
	z_b	(cm)					#	
	z_{MI}	(cm)	4.50					
	$z_{pii,\alpha}$	(cm)	4.50					
	f_{awf}	(MHz)	2.47	3.46		3.46		3.03
Other Information	p_{rr}	(Hz)	#					
	s_{rr}	(Hz)	77.37					
	n_{pps}	-	8.00					
	$I_{pa,\alpha}$ at $z_{pii,\alpha}$	(W/cm ²)	448.17					
	$I_{spta,\alpha}$ at $z_{pii,\alpha}$ or $z_{sii,\alpha}$	(mW/cm ²)	32.62					
	I_{spta} at z_{pii} or z_{sii}	(mW/cm ²)	67.94					
	p_r at z_{pii}	(MPa)	3.37					
Operating Control Conditions	TX Level	(dB)	0.00	0.00		0.00		0.00
	Focus	(mm)	60.00	40.00		40.00		20.00
a This Index is not relevant to this operating mode.								
b This transducer is not intended for transcranial or neonatal cephalic uses.								
c This formulation for TIS is less than that for an alternate formulation in this mode.								
# No data is provided for this operation condition since the maximum index value is not reported for the reason listed.								

Acoustic Output Reporting Table — Track 3, FDA 510(k) and IEC 60601-2-37

(Per transducer/mode that exceeds MI or TI value of 1.0)

Transducer model: **5P1**Operating mode: **Steerable Continuous Wave Doppler**

Index Label			MI	TIS		TIB		TIC
				At surface	Below surface	At surface	Below surface	
Maximum Index Value			0.09	1.11		3.47		3.10
Index Component Value				1.02	1.11	1.02	3.47	
Acoustic Parameters	IEC	Units						
	$p_{ra,\alpha}$ at z_{MI}	(MPa)	0.12					
	P	(mW)		#		165.23		165.23
	$P_{1\times1}$	(mW)		118.25		#		
	z_s	(cm)			2.00			
	z_b	(cm)					3.50	
	z_{MI}	(cm)	3.20					
	$z_{pII,\alpha}$	(cm)	3.20					
	f_{awf}	(MHz)	1.80	1.80		1.80		1.80
Other Information	p_{rr}	(Hz)	1818200.00					
	s_{rr}	(Hz)	#					
	n_{pps}	-	#					
	$I_{pa,\alpha}$ at $z_{pII,\alpha}$	(W/cm ²)	0.50					
	$I_{spta,\alpha}$ at $z_{pII,\alpha}$ or $z_{sII,\alpha}$	(mW/cm ²)	502.36					
	I_{spta} at z_{pII} or z_{sII}	(mW/cm ²)	766.55					
	p_r at z_{pII}	(MPa)	0.15					
Operating Control Conditions	TX Level	(dB)	0.00	0.00		0.00		0.00
	Focus	(mm)	40.00	120.00		120.00		120.00
a	This Index is not relevant to this operating mode.							
b	This transducer is not intended for transcranial or neonatal cephalic uses.							
c	This formulation for TIS is less than that for an alternate formulation in this mode.							
#	No data is provided for this operation condition since the maximum index value is not reported for the reason listed.							

Acoustic Output Reporting Table — Track 3, FDA 510(k) and IEC 60601-2-37

(Per transducer/mode that exceeds MI or TI value of 1.0)

Transducer model: **8V4**Operating mode: **2D-mode**

Index Label			MI	TIS		TIB		TIC
				At surface	Below surface	At surface	Below surface	
Maximum Index Value			1.50	0.97		0.97		1.10
Index Component Value				0.97	0.97	0.97	0.97	
Acoustic Parameters	IEC	Units						
	$p_{ra,\alpha}$ at z_{MI}	(MPa)	3.41					
	P	(mW)		#		36.19		43.55
	P_{1x1}	(mW)		36.19		36.19		
	z_s	(cm)			#			
	z_b	(cm)					#	
	z_{MI}	(cm)	0.50					
	$z_{pii,\alpha}$	(cm)	1.50					
	f_{awf}	(MHz)	4.63	5.66		5.66		3.91
Other Information	p_{rr}	(Hz)	#					
	s_{rr}	(Hz)	30.00					
	n_{pps}	-	1.00					
	$I_{pa,\alpha}$ at $z_{pii,\alpha}$	(W/cm ²)	183.10					
	$I_{spta,\alpha}$ at $z_{pii,\alpha}$ or $z_{sii,\alpha}$	(mW/cm ²)	50.30					
	I_{spta} at z_{pii} or z_{sii}	(mW/cm ²)	81.19					
	p_r at z_{pii}	(MPa)	3.11					
Operating Control Conditions	TX Level	(dB)	0.00	0.00		0.00		0.00
	Focus	(mm)	80.00	20.00		20.00		30.00
a This Index is not relevant to this operating mode.								
b This transducer is not intended for transcranial or neonatal cephalic uses.								
c This formulation for TIS is less than that for an alternate formulation in this mode.								
# No data is provided for this operation condition since the maximum index value is not reported for the reason listed.								

Acoustic Output Reporting Table — Track 3, FDA 510(k) and IEC 60601-2-37

(Per transducer/mode that exceeds MI or TI value of 1.0)

Transducer model: **8V4**Operating mode: **M-mode**

Index Label			MI	TIS		TIB		TIC
				At surface	Below surface	At surface	Below surface	
Maximum Index Value			1.49	0.82		1.15		1.00
Index Component Value				0.82	0.50	0.73	1.15	
Acoustic Parameters	IEC	Units						
	$p_{ra,\alpha}$ at z_{MI}	(MPa)	2.68					
	P	(mW)		#		39.53		39.53
	$P_{1\times1}$	(mW)		36.27		39.53		
	z_s	(cm)			1.50			
	z_b	(cm)					2.70	
	z_{MI}	(cm)	2.70					
	$z_{pII,\alpha}$	(cm)	2.70					
	f_{awf}	(MHz)	3.29	4.73		3.90		3.90
Other Information	p_{rr}	(Hz)	720.00					
	s_{rr}	(Hz)	#					
	n_{pps}	-	#					
	$I_{pa,\alpha}$ at $z_{pII,\alpha}$	(W/cm ²)	157.05					
	$I_{spta,\alpha}$ at $z_{pII,\alpha}$ or $z_{sII,\alpha}$	(mW/cm ²)	90.33					
	I_{spta} at z_{pII} or z_{sII}	(mW/cm ²)	168.19					
	p_r at z_{pII}	(MPa)	3.44					
Operating Control Conditions	TX Level	(dB)	0.00	0.00		0.00		0.00
	Focus	(mm)	30.00	30.00		40.00		40.00
a	This Index is not relevant to this operating mode.							
b	This transducer is not intended for transcranial or neonatal cephalic uses.							
c	This formulation for TIS is less than that for an alternate formulation in this mode.							
#	No data is provided for this operation condition since the maximum index value is not reported for the reason listed.							

Acoustic Output Reporting Table — Track 3, FDA 510(k) and IEC 60601-2-37

(Per transducer/mode that exceeds MI or TI value of 1.0)

Transducer model: **8V4**Operating mode: **Pulsed Wave Doppler**

Index Label			MI	TIS		TIB		TIC
				At surface	Below surface	At surface	Below surface	
Maximum Index Value			1.41	0.85		1.93		0.99
Index Component Value				0.85	0.54	0.52	1.93	
Acoustic Parameters	IEC	Units						
	$p_{ra,\alpha}$ at z_{MI}	(MPa)	2.95					
	P	(mW)		#		30.45		39.35
	P_{1x1}	(mW)		41.66		30.45		
	z_s	(cm)			1.50			
	z_b	(cm)					1.80	
	z_{MI}	(cm)	1.00					
	$z_{pii,\alpha}$	(cm)	1.50					
	f_{awf}	(MHz)	4.23	4.28		3.62		4.24
Other Information	p_{rr}	(Hz)	243.00					
	s_{rr}	(Hz)	#					
	n_{pps}	-	#					
	$I_{pa,\alpha}$ at $z_{pii,\alpha}$	(W/cm ²)	201.77					
	$I_{spta,\alpha}$ at $z_{pii,\alpha}$ or $z_{sii,\alpha}$	(mW/cm ²)	42.58					
	I_{spta} at z_{pii} or z_{sii}	(mW/cm ²)	107.19					
	p_r at z_{pii}	(MPa)	3.17					
Operating Control Conditions	TX Level	(dB)	0.00	0.00		0.00		0.00
	Focus	(mm)	60.00	20.00		20.00		30.00
a This Index is not relevant to this operating mode.								
b This transducer is not intended for transcranial or neonatal cephalic uses.								
c This formulation for TIS is less than that for an alternate formulation in this mode.								
# No data is provided for this operation condition since the maximum index value is not reported for the reason listed.								

Acoustic Output Reporting Table — Track 3, FDA 510(k) and IEC 60601-2-37

(Per transducer/mode that exceeds MI or TI value of 1.0)

Transducer model: **8V4**Operating mode: **Color Flow Doppler**

Index Label			MI	TIS		TIB		TIC
				At surface	Below surface	At surface	Below surface	
Maximum Index Value			1.49	0.89		0.89		1.06
Index Component Value				0.41	0.41	0.41	0.41	
Acoustic Parameters	IEC	Units						
	$p_{ra,\alpha}$ at z_{MI}	(MPa)	2.85					
	P	(mW)		#		20.03		20.04
	P_{1x1}	(mW)		20.03		20.03		
	z_s	(cm)			#			
	z_b	(cm)					#	
	z_{MI}	(cm)	2.40					
	$z_{pii,\alpha}$	(cm)	2.40					
	f_{awf}	(MHz)	3.62	4.28		4.28		4.22
Other Information	pr_r	(Hz)	#					
	sr_r	(Hz)	77.84					
	n_{pps}	-	10.00					
	$I_{pa,\alpha}$ at $z_{pii,\alpha}$	(W/cm ²)	427.14					
	$I_{spta,\alpha}$ at $z_{pii,\alpha}$ or $z_{sii,\alpha}$	(mW/cm ²)	3.85					
	I_{spta} at z_{pii} or z_{sii}	(mW/cm ²)	13.78					
	p_r at z_{pii}	(MPa)	3.84					
Operating Control Conditions	TX Level	(dB)	0.00	0.00		0.00		0.00
	Focus	(mm)	30.00	20.00		20.00		100.00
a This Index is not relevant to this operating mode.								
b This transducer is not intended for transcranial or neonatal cephalic uses.								
c This formulation for TIS is less than that for an alternate formulation in this mode.								
# No data is provided for this operation condition since the maximum index value is not reported for the reason listed.								

Acoustic Output Reporting Table — Track 3, FDA 510(k) and IEC 60601-2-37

(Per transducer/mode that exceeds MI or TI value of 1.0)

Transducer model: **8V4**Operating mode: **Steerable Continuous Wave Doppler**

Index Label			MI	TIS		TIB		TIC
				At surface	Below surface	At surface	Below surface	
Maximum Index Value			0.06	0.46		1.71		0.91
Index Component Value				0.46	0.35	0.46	1.71	
Acoustic Parameters	IEC	Units						
	$p_{ra,\alpha}$ at z_{MI}	(MPa)	0.11					
	P	(mW)		#		24.31		24.31
	P_{1x1}	(mW)		24.31		24.31		
	z_s	(cm)			1.00			
	z_b	(cm)					1.20	
	z_{MI}	(cm)	1.30					
	$z_{pii,\alpha}$	(cm)	1.30					
	f_{awf}	(MHz)	4.00	4.00		4.00		4.00
Other Information	p_{rr}	(Hz)	4000000.00					
	s_{rr}	(Hz)	#					
	n_{pps}	-	#					
	$I_{pa,\alpha}$ at $z_{pii,\alpha}$	(W/cm ²)	0.42					
	$I_{spta,\alpha}$ at $z_{pii,\alpha}$ or $z_{sii,\alpha}$	(mW/cm ²)	416.84					
	I_{spta} at z_{pii} or z_{sii}	(mW/cm ²)	646.31					
	p_r at z_{pii}	(MPa)	0.14					
Operating Control Conditions	TX Level	(dB)	0.00	0.00		0.00		0.00
	Focus	(mm)	20.00	20.00		20.00		20.00
a This Index is not relevant to this operating mode.								
b This transducer is not intended for transcranial or neonatal cephalic uses.								
c This formulation for TIS is less than that for an alternate formulation in this mode.								
# No data is provided for this operation condition since the maximum index value is not reported for the reason listed.								

Acoustic Output Reporting Table — Track 3, FDA 510(k) and IEC 60601-2-37

(Per transducer/mode that exceeds MI or TI value of 1.0)

Transducer model: **10V4**Operating mode: **2D-mode**

Index Label			MI	TIS		TIB		TIC
				At surface	Below surface	At surface	Below surface	
Maximum Index Value			1.47	0.88		0.88		1.77
Index Component Value				0.88	0.88	0.88	0.88	
Acoustic Parameters	IEC	Units						
	$p_{ra,\alpha}$ at z_{MI}	(MPa)	3.40					
	P	(mW)		#		49.81		30.88
	$P_{1\times1}$	(mW)		45.97		45.97		
	z_s	(cm)			#			
	z_b	(cm)					#	
	z_{MI}	(cm)	0.70					
	$z_{pii,\alpha}$	(cm)	1.80					
	f_{awf}	(MHz)	5.20	4.03		4.03		4.04
Other Information	p_{rr}	(Hz)	#					
	s_{rr}	(Hz)	79.97					
	n_{pps}	-	1.00					
	$I_{pa,\alpha}$ at $z_{pii,\alpha}$	(W/cm ²)	203.08					
	$I_{spta,\alpha}$ at $z_{pii,\alpha}$ or $z_{sii,\alpha}$	(mW/cm ²)	107.86					
	I_{spta} at z_{pii} or z_{sii}	(mW/cm ²)	198.46					
	p_r at z_{pii}	(MPa)	3.58					
	Operating Control Conditions	TX Level	(dB)	0.00	0.00		0.00	
Focus		(mm)	65.00	30.00		30.00		5.00
a	This Index is not relevant to this operating mode.							
b	This transducer is not intended for transcranial or neonatal cephalic uses.							
c	This formulation for TIS is less than that for an alternate formulation in this mode.							
#	No data is provided for this operation condition since the maximum index value is not reported for the reason listed.							

Acoustic Output Reporting Table — Track 3, FDA 510(k) and IEC 60601-2-37

(Per transducer/mode that exceeds MI or TI value of 1.0)

Transducer model: **10V4**Operating mode: **M-mode**

Index Label			MI	TIS		TIB		TIC
				At surface	Below surface	At surface	Below surface	
Maximum Index Value			1.47	0.69		0.78		0.70
Index Component Value				0.69	0.33	0.41	0.78	
Acoustic Parameters	IEC	Units						
	$p_{ra,\alpha}$ at z_{MI}	(MPa)	3.40					
	P	(mW)		#		21.44		27.33
	P_{1x1}	(mW)		20.54		21.44		
	z_s	(cm)			1.50			
	z_b	(cm)					2.50	
	z_{MI}	(cm)	0.70					
	$z_{pii,\alpha}$	(cm)	1.80					
	f_{awf}	(MHz)	5.20	7.08		4.03		5.17
Other Information	p_{rr}	(Hz)	720.00					
	s_{rr}	(Hz)	#					
	n_{pps}	-	#					
	$I_{pa,\alpha}$ at $z_{pii,\alpha}$	(W/cm ²)	203.08					
	$I_{spta,\alpha}$ at $z_{pii,\alpha}$ or $z_{sii,\alpha}$	(mW/cm ²)	64.86					
	I_{spta} at z_{pii} or z_{sii}	(mW/cm ²)	135.76					
	p_r at z_{pii}	(MPa)	3.58					
Operating Control Conditions	TX Level	(dB)	0.00	0.00		0.00		0.00
	Focus	(mm)	65.00	65.00		30.00		30.00
a This Index is not relevant to this operating mode.								
b This transducer is not intended for transcranial or neonatal cephalic uses.								
c This formulation for TIS is less than that for an alternate formulation in this mode.								
# No data is provided for this operation condition since the maximum index value is not reported for the reason listed.								

Acoustic Output Reporting Table — Track 3, FDA 510(k) and IEC 60601-2-37

(Per transducer/mode that exceeds MI or TI value of 1.0)

Transducer model: **10V4**Operating mode: **Pulsed Wave Doppler**

Index Label			MI	TIS		TIB		TIC
				At surface	Below surface	At surface	Below surface	
Maximum Index Value			1.52	1.38		1.80		1.38
Index Component Value				1.38	0.69	0.46	1.80	
Acoustic Parameters	IEC	Units						
	$p_{ra,\alpha}$ at z_{MI}	(MPa)	3.41					
	P	(mW)		#		23.78		28.12
	$P_{1\times1}$	(mW)		48.95		23.78		
	z_s	(cm)			1.80			
	z_b	(cm)					1.10	
	z_{MI}	(cm)	0.70					
	$z_{pii,\alpha}$	(cm)	2.10					
	f_{awf}	(MHz)	5.02	5.93		4.06		6.02
Other Information	p_{rr}	(Hz)	243.00					
	s_{rr}	(Hz)	#					
	n_{pps}	-	#					
	$I_{pa,\alpha}$ at $z_{pii,\alpha}$	(W/cm ²)	205.65					
	$I_{spta,\alpha}$ at $z_{pii,\alpha}$ or $z_{sii,\alpha}$	(mW/cm ²)	33.94					
	I_{spta} at z_{pii} or z_{sii}	(mW/cm ²)	67.41					
	p_r at z_{pii}	(MPa)	3.27					
Operating Control Conditions	TX Level	(dB)	0.00	0.00		0.00		0.00
	Focus	(mm)	90.00	65.00		10.00		5.00
a	This Index is not relevant to this operating mode.							
b	This transducer is not intended for transcranial or neonatal cephalic uses.							
c	This formulation for TIS is less than that for an alternate formulation in this mode.							
#	No data is provided for this operation condition since the maximum index value is not reported for the reason listed.							

Acoustic Output Reporting Table — Track 3, FDA 510(k) and IEC 60601-2-37

(Per transducer/mode that exceeds MI or TI value of 1.0)

Transducer model: **10V4**Operating mode: **Color Flow Doppler**

Index Label			MI	TIS		TIB		TIC
				At surface	Below surface	At surface	Below surface	
Maximum Index Value			1.49	1.20		1.20		1.50
Index Component Value				1.15	1.15	1.15	1.15	
Acoustic Parameters	IEC	Units						
	$p_{ra,\alpha}$ at z_{MI}	(MPa)	3.35					
	P	(mW)		#		41.44		18.61
	P_{1x1}	(mW)		41.34		41.34		
	z_s	(cm)			#			
	z_b	(cm)					#	
	z_{MI}	(cm)	1.90					
	$z_{pii,\alpha}$	(cm)	1.90					
	f_{awf}	(MHz)	4.96	5.86		5.86		6.02
Other Information	p_{rr}	(Hz)	#					
	s_{rr}	(Hz)	353.97					
	n_{pps}	-	11.00					
	$I_{pa,\alpha}$ at $z_{pii,\alpha}$	(W/cm ²)	657.08					
	$I_{spta,\alpha}$ at $z_{pii,\alpha}$ or $z_{sii,\alpha}$	(mW/cm ²)	128.99					
	I_{spta} at z_{pii} or z_{sii}	(mW/cm ²)	223.13					
	p_r at z_{pii}	(MPa)	4.64					
Operating Control Conditions	TX Level	(dB)	0.00	0.00		0.00		0.00
	Focus	(mm)	20.00	20.00		20.00		5.00
a This Index is not relevant to this operating mode.								
b This transducer is not intended for transcranial or neonatal cephalic uses.								
c This formulation for TIS is less than that for an alternate formulation in this mode.								
# No data is provided for this operation condition since the maximum index value is not reported for the reason listed.								

Acoustic Output Reporting Table — Track 3, FDA 510(k) and IEC 60601-2-37

(Per transducer/mode that exceeds MI or TI value of 1.0)

Transducer model: **10V4**Operating mode: **Steerable Continuous Wave Doppler**

Index Label			MI	TIS		TIB		TIC
				At surface	Below surface	At surface	Below surface	
Maximum Index Value			0.06	1.05		2.12		1.57
Index Component Value				1.05	0.73	0.67	2.12	
Acoustic Parameters	IEC	Units						
	$p_{ra,\alpha}$ at z_{MI}	(MPa)	0.12					
	P	(mW)		#		31.78		49.56
	$P_{1\times1}$	(mW)		49.56		31.78		
	z_s	(cm)			1.20			
	z_b	(cm)					1.10	
	z_{MI}	(cm)	1.30					
	$z_{pii,\alpha}$	(cm)	1.30					
	f_{awf}	(MHz)	4.00	4.44		4.44		4.44
Other Information	p_{rr}	(Hz)	4000000.00					
	s_{rr}	(Hz)	#					
	n_{pps}	-	#					
	$I_{pa,\alpha}$ at $z_{pii,\alpha}$	(W/cm ²)	0.50					
	$I_{spta,\alpha}$ at $z_{pii,\alpha}$ or $z_{sii,\alpha}$	(mW/cm ²)	500.50					
	I_{spta} at z_{pii} or z_{sii}	(mW/cm ²)	757.21					
	p_r at z_{pii}	(MPa)	0.15					
Operating Control Conditions	TX Level	(dB)	0.00	0.00		0.00		0.00
	Focus	(mm)	15.00	90.00		10.00		90.00
a	This Index is not relevant to this operating mode.							
b	This transducer is not intended for transcranial or neonatal cephalic uses.							
c	This formulation for TIS is less than that for an alternate formulation in this mode.							
#	No data is provided for this operation condition since the maximum index value is not reported for the reason listed.							

Acoustic Output Reporting Table — Track 3, FDA 510(k) and IEC 60601-2-37

(Per transducer/mode that exceeds MI or TI value of 1.0)

Transducer model: **5VT**Operating mode: **2D-mode**

Index Label			MI	TIS		TIB		TIC
				At surface	Below surface	At surface	Below surface	
Maximum Index Value			1.14	0.61		0.61		0.83
Index Component Value				0.61	0.61	0.61	0.61	
Acoustic Parameters	IEC	Units						
	$p_{ra,\alpha}$ at z_{MI}	(MPa)	2.30					
	P	(mW)		#		24.45		28.40
	P_{1x1}	(mW)		24.45		24.45		
	z_s	(cm)			#			
	z_b	(cm)					#	
	z_{MI}	(cm)	1.90					
	$z_{pii,\alpha}$	(cm)	1.90					
	f_{awf}	(MHz)	4.02	5.28		5.28		4.22
Other Information	p_{rr}	(Hz)	#					
	s_{rr}	(Hz)	149.92					
	n_{pps}	-	1.00					
	$I_{pa,\alpha}$ at $z_{pii,\alpha}$	(W/cm ²)	256.85					
	$I_{spta,\alpha}$ at $z_{pii,\alpha}$ or $z_{sii,\alpha}$	(mW/cm ²)	156.13					
	I_{spta} at z_{pii} or z_{sii}	(mW/cm ²)	243.24					
	p_r at z_{pii}	(MPa)	2.99					
Operating Control Conditions	TX Level	(dB)	0.00	0.00		0.00		0.00
	Focus	(mm)	20.00	20.00		20.00		20.00
a	This Index is not relevant to this operating mode.							
b	This transducer is not intended for transcranial or neonatal cephalic uses.							
c	This formulation for TIS is less than that for an alternate formulation in this mode.							
#	No data is provided for this operation condition since the maximum index value is not reported for the reason listed.							

Acoustic Output Reporting Table — Track 3, FDA 510(k) and IEC 60601-2-37

(Per transducer/mode that exceeds MI or TI value of 1.0)

Transducer model: **5VT**Operating mode: **M-mode**

Index Label			MI	TIS		TIB		TIC
				At surface	Below surface	At surface	Below surface	
Maximum Index Value			1.45	0.59		1.07		0.79
Index Component Value				0.59	0.39	0.51	1.07	
Acoustic Parameters	IEC	Units						
	$p_{ra,\alpha}$ at z_{MI}	(MPa)	2.92					
	P	(mW)		#		25.27		26.73
	P_{1x1}	(mW)		25.97		25.27		
	z_s	(cm)			1.30			
	z_b	(cm)					1.80	
	z_{MI}	(cm)	1.90					
	$z_{pii,\alpha}$	(cm)	1.90					
	f_{awf}	(MHz)	4.02	4.80		4.23		4.22
Other Information	prr	(Hz)	720.00					
	srr	(Hz)	#					
	n_{pps}	-	#					
	$I_{pa,\alpha}$ at $z_{pii,\alpha}$	(W/cm ²)	444.92					
	$I_{spta,\alpha}$ at $z_{pii,\alpha}$ or $z_{sii,\alpha}$	(mW/cm ²)	201.12					
	I_{spta} at z_{pii} or z_{sii}	(mW/cm ²)	340.49					
	p_r at z_{pii}	(MPa)	3.80					
Operating Control Conditions	TX Level	(dB)	0.00	0.00		0.00		0.00
	Focus	(mm)	20.00	20.00		20.00		30.00
a This Index is not relevant to this operating mode.								
b This transducer is not intended for transcranial or neonatal cephalic uses.								
c This formulation for TIS is less than that for an alternate formulation in this mode.								
# No data is provided for this operation condition since the maximum index value is not reported for the reason listed.								

Acoustic Output Reporting Table — Track 3, FDA 510(k) and IEC 60601-2-37

(Per transducer/mode that exceeds MI or TI value of 1.0)

Transducer model: **5VT**Operating mode: **Pulsed Wave Doppler**

Index Label			MI	TIS		TIB		TIC
				At surface	Below surface	At surface	Below surface	
Maximum Index Value			1.46	0.46		1.37		0.74
Index Component Value				0.46	0.33	0.48	1.37	
Acoustic Parameters	IEC	Units						
	$p_{ra,\alpha}$ at z_{MI}	(MPa)	2.59					
	P	(mW)		#		26.52		25.28
	P_{1x1}	(mW)		25.28		26.52		
	z_s	(cm)			1.30			
	z_b	(cm)					1.70	
	z_{MI}	(cm)	2.50					
	$z_{pii,\alpha}$	(cm)	2.50					
	f_{awf}	(MHz)	3.19	3.83		3.84		3.83
Other Information	p_{rr}	(Hz)	243.00					
	s_{rr}	(Hz)	#					
	n_{pps}	-	#					
	$I_{pa,\alpha}$ at $z_{pii,\alpha}$	(W/cm ²)	405.15					
	$I_{spta,\alpha}$ at $z_{pii,\alpha}$ or $z_{sii,\alpha}$	(mW/cm ²)	128.32					
	I_{spta} at z_{pii} or z_{sii}	(mW/cm ²)	217.73					
	p_r at z_{pii}	(MPa)	3.29					
Operating Control Conditions	TX Level	(dB)	0.00	0.00		0.00		0.00
	Focus	(mm)	30.00	50.00		20.00		50.00
a	This Index is not relevant to this operating mode.							
b	This transducer is not intended for transcranial or neonatal cephalic uses.							
c	This formulation for TIS is less than that for an alternate formulation in this mode.							
#	No data is provided for this operation condition since the maximum index value is not reported for the reason listed.							

Acoustic Output Reporting Table — Track 3, FDA 510(k) and IEC 60601-2-37

(Per transducer/mode that exceeds MI or TI value of 1.0)

Transducer model: **5VT**Operating mode: **Color Flow Doppler**

Index Label			MI	TIS		TIB		TIC
				At surface	Below surface	At surface	Below surface	
Maximum Index Value			1.46	0.24		0.24		0.39
Index Component Value				0.24	0.24	0.24	0.24	
Acoustic Parameters	IEC	Units						
	$p_{ra,\alpha}$ at z_{MI}	(MPa)	2.39					
	P	(mW)		#		13.22		13.22
	P_{1x1}	(mW)		13.22		13.22		
	z_s	(cm)			#			
	z_b	(cm)					#	
	z_{MI}	(cm)	2.00					
	$z_{pii,\alpha}$	(cm)	2.00					
	f_{awf}	(MHz)	2.69	3.84		3.84		3.81
Other Information	pr_r	(Hz)	#					
	sr_r	(Hz)	15.00					
	n_{pps}	-	2.00					
	$I_{pa,\alpha}$ at $z_{pii,\alpha}$	(W/cm ²)	124.91					
	$I_{spta,\alpha}$ at $z_{pii,\alpha}$ or $z_{sii,\alpha}$	(mW/cm ²)	7.04					
	I_{spta} at z_{pii} or z_{sii}	(mW/cm ²)	9.65					
	p_r at z_{pii}	(MPa)	2.80					
Operating Control Conditions	TX Level	(dB)	0.00	0.00		0.00		0.00
	Focus	(mm)	20.00	20.00		20.00		160.00
a This Index is not relevant to this operating mode.								
b This transducer is not intended for transcranial or neonatal cephalic uses.								
c This formulation for TIS is less than that for an alternate formulation in this mode.								
# No data is provided for this operation condition since the maximum index value is not reported for the reason listed.								

Acoustic Output Reporting Table — Track 3, FDA 510(k) and IEC 60601-2-37

(Per transducer/mode that exceeds MI or TI value of 1.0)

Transducer model: **P4-2**Operating mode: **2D-mode**

Index Label			MI	TIS		TIB		TIC
				At surface	Below surface	At surface	Below surface	
Maximum Index Value			1.41	1.45		1.45		2.93
Index Component Value				1.45	1.45	1.45	1.45	
Acoustic Parameters	IEC	Units						
	$p_{ra,\alpha}$ at z_{MI}	(MPa)	1.90					
	P	(mW)		#		130.84		133.48
	P_{1x1}	(mW)		96.92		96.92		
	z_s	(cm)			#			
	z_b	(cm)					#	
	z_{MI}	(cm)	1.80					
	$z_{pii,\alpha}$	(cm)	1.80					
	f_{awf}	(MHz)	1.82	3.15		3.15		2.03
Other Information	p_{rr}	(Hz)	#					
	s_{rr}	(Hz)	46.13					
	η_{pps}	-	2.00					
	$I_{pa,\alpha}$ at $z_{pii,\alpha}$	(W/cm ²)	104.54					
	$I_{spta,\alpha}$ at $z_{pii,\alpha}$ or $z_{sii,\alpha}$	(mW/cm ²)	12.86					
	I_{spta} at z_{pii} or z_{sii}	(mW/cm ²)	23.78					
	p_r at z_{pii}	(MPa)	2.13					
Operating Control Conditions	TX Level	(dB)	0.00	0.00		0.00		0.00
	Focus	(mm)	30.00	50.00		50.00		30.00
a This Index is not relevant to this operating mode.								
b This transducer is not intended for transcranial or neonatal cephalic uses.								
c This formulation for TIS is less than that for an alternate formulation in this mode.								
# No data is provided for this operation condition since the maximum index value is not reported for the reason listed.								

Acoustic Output Reporting Table — Track 3, FDA 510(k) and IEC 60601-2-37

(Per transducer/mode that exceeds MI or TI value of 1.0)

Transducer model: **P4-2**Operating mode: **M-mode**

Index Label			MI	TIS		TIB		TIC
				At surface	Below surface	At surface	Below surface	
Maximum Index Value			1.50	0.93		3.78		2.64
Index Component Value				0.93	0.84	0.89	3.78	
Acoustic Parameters	IEC	Units						
	$p_{ra,\alpha}$ at z_{MI}	(MPa)	2.43					
	P	(mW)		#		148.43		137.87
	$P_{1\times1}$	(mW)		73.61		93.97		
	z_s	(cm)			4.90			
	z_b	(cm)					4.10	
	z_{MI}	(cm)	2.80					
	$z_{pii,\alpha}$	(cm)	2.80					
	f_{awf}	(MHz)	2.51	2.66		2.00		2.00
Other Information	prr	(Hz)	720.00					
	srr	(Hz)	#					
	n_{pps}	-	#					
	$I_{pa,\alpha}$ at $z_{pii,\alpha}$	(W/cm ²)	269.52					
	$I_{spta,\alpha}$ at $z_{pii,\alpha}$ or $z_{sii,\alpha}$	(mW/cm ²)	177.93					
	I_{spta} at z_{pii} or z_{sii}	(mW/cm ²)	288.62					
	p_r at z_{pii}	(MPa)	2.90					
	Operating Control Conditions	TX Level	(dB)	0.00	0.00		0.00	
Focus		(mm)	50.00	110.00		70.00		50.00
a	This Index is not relevant to this operating mode.							
b	This transducer is not intended for transcranial or neonatal cephalic uses.							
c	This formulation for TIS is less than that for an alternate formulation in this mode.							
#	No data is provided for this operation condition since the maximum index value is not reported for the reason listed.							

Acoustic Output Reporting Table — Track 3, FDA 510(k) and IEC 60601-2-37

(Per transducer/mode that exceeds MI or TI value of 1.0)

Transducer model: **P4-2**Operating mode: **Pulsed Wave Doppler**

Index Label			MI	TIS		TIB		TIC
				At surface	Below surface	At surface	Below surface	
Maximum Index Value			1.50	1.65		4.25		2.78
Index Component Value				1.65	1.12	1.16	4.25	
Acoustic Parameters	IEC	Units						
	$p_{ra,\alpha}$ at z_{MI}	(MPa)	2.23					
	P	(mW)		#		122.82		188.18
	P_{1x1}	(mW)		107.61		120.66		
	z_s	(cm)			1.80			
	z_b	(cm)					1.80	
	z_{MI}	(cm)	5.10					
	$z_{pii,\alpha}$	(cm)	5.10					
	f_{awf}	(MHz)	2.23	3.21		2.02		2.02
Other Information	p_{rr}	(Hz)	610.35					
	s_{rr}	(Hz)	#					
	n_{pps}	-	#					
	$I_{pa,\alpha}$ at $z_{pii,\alpha}$	(W/cm ²)	372.57					
	$I_{spta,\alpha}$ at $z_{pii,\alpha}$ or $z_{sii,\alpha}$	(mW/cm ²)	304.81					
	I_{spta} at z_{pii} or z_{sii}	(mW/cm ²)	664.96					
	p_r at z_{pii}	(MPa)	3.15					
Operating Control Conditions	TX Level	(dB)	0.00	0.00		0.00		0.00
	Focus	(mm)	70.00	30.00		30.00		110.00
a This Index is not relevant to this operating mode.								
b This transducer is not intended for transcranial or neonatal cephalic uses.								
c This formulation for TIS is less than that for an alternate formulation in this mode.								
# No data is provided for this operation condition since the maximum index value is not reported for the reason listed.								

Acoustic Output Reporting Table — Track 3, FDA 510(k) and IEC 60601-2-37

(Per transducer/mode that exceeds MI or TI value of 1.0)

Transducer model: **P4-2**Operating mode: **Color Flow Doppler**

Index Label			MI	TIS		TIB		TIC
				At surface	Below surface	At surface	Below surface	
Maximum Index Value			1.50	1.34		1.34		2.86
Index Component Value				0.62	0.62	0.62	0.62	
Acoustic Parameters	IEC	Units						
	$p_{ra,\alpha}$ at z_{MI}	(MPa)	2.02					
	P	(mW)		#		54.86		67.69
	P_{1x1}	(mW)		40.64		40.64		
	z_s	(cm)			#			
	z_b	(cm)					#	
	z_{MI}	(cm)	5.60					
	$z_{pII,\alpha}$	(cm)	5.60					
	f_{awf}	(MHz)	1.83	3.21		3.21		2.02
Other Information	p_{rr}	(Hz)	#					
	s_{rr}	(Hz)	14.48					
	n_{pps}	-	2.00					
	$I_{pa,\alpha}$ at $z_{pII,\alpha}$	(W/cm ²)	208.80					
	$I_{spta,\alpha}$ at $z_{pII,\alpha}$ or $z_{sII,\alpha}$	(mW/cm ²)	7.40					
	I_{spta} at z_{pII} or z_{sII}	(mW/cm ²)	13.73					
	p_r at z_{pII}	(MPa)	2.76					
Operating Control Conditions	TX Level	(dB)	0.00	0.00		0.00		0.00
	Focus	(mm)	90.00	30.00		30.00		30.00
a	This Index is not relevant to this operating mode.							
b	This transducer is not intended for transcranial or neonatal cephalic uses.							
c	This formulation for TIS is less than that for an alternate formulation in this mode.							
#	No data is provided for this operation condition since the maximum index value is not reported for the reason listed.							

Acoustic Output Reporting Table — Track 3, FDA 510(k) and IEC 60601-2-37

(Per transducer/mode that exceeds MI or TI value of 1.0)

Transducer model: **P4-2**Operating mode: **Steerable Continuous Wave Doppler**

Index Label			MI	TIS		TIB		TIC
				At surface	Below surface	At surface	Below surface	
Maximum Index Value			0.09	0.76		3.82		2.06
Index Component Value				0.76	0.68	0.75	3.82	
Acoustic Parameters	IEC	Units						
	$p_{ra,\alpha}$ at z_{MI}	(MPa)	0.12					
	P	(mW)		#		96.94		98.32
	$P_{1\times1}$	(mW)		87.54		86.30		
	z_s	(cm)			1.80			
	z_b	(cm)					1.90	
	z_{MI}	(cm)	2.00					
	$z_{pII,\alpha}$	(cm)	2.00					
	f_{awf}	(MHz)	1.82	1.82		1.82		1.82
Other Information	p_{rr}	(Hz)	1818200.00					
	s_{rr}	(Hz)	#					
	n_{pps}	-	#					
	$I_{pa,\alpha}$ at $z_{pII,\alpha}$	(W/cm ²)	0.45					
	$I_{spta,\alpha}$ at $z_{pII,\alpha}$ or $z_{sII,\alpha}$	(mW/cm ²)	450.21					
	I_{spta} at z_{pII} or z_{sII}	(mW/cm ²)	668.95					
	p_r at z_{pII}	(MPa)	0.14					
Operating Control Conditions	TX Level	(dB)	0.00	0.00		0.00		0.00
	Focus	(mm)	110.00	180.00		150.00		180.00
a	This Index is not relevant to this operating mode.							
b	This transducer is not intended for transcranial or neonatal cephalic uses.							
c	This formulation for TIS is less than that for an alternate formulation in this mode.							
#	No data is provided for this operation condition since the maximum index value is not reported for the reason listed.							

Acoustic Output Reporting Table — Track 3, FDA 510(k) and IEC 60601-2-37

(Per transducer/mode that exceeds MI or TI value of 1.0)

Transducer model: **9MC3**Operating mode: **2D-mode**

Index Label			MI	TIS		TIB		TIC
				At surface	Below surface	At surface	Below surface	
Maximum Index Value			1.50	0.95		0.95		0.82
Index Component Value				0.95	0.95	0.95	0.95	
Acoustic Parameters	IEC	Units						
	$p_{ra,\alpha}$ at z_{MI}	(MPa)	3.21					
	P	(mW)		#		38.15		45.29
	$P_{1\times1}$	(mW)		34.28		34.28		
	z_s	(cm)			#			
	z_b	(cm)					#	
	z_{MI}	(cm)	1.40					
	$z_{pII,\alpha}$	(cm)	1.40					
	f_{awf}	(MHz)	4.56	5.81		5.81		4.57
Other Information	pr_r	(Hz)	#					
	sr_r	(Hz)	21.32					
	n_{pps}	-	1.00					
	$I_{pa,\alpha}$ at $z_{pII,\alpha}$	(W/cm ²)	603.80					
	$I_{spta,\alpha}$ at $z_{pII,\alpha}$ or $z_{sII,\alpha}$	(mW/cm ²)	19.39					
	I_{spta} at z_{pII} or z_{sII}	(mW/cm ²)	30.11					
	p_r at z_{pII}	(MPa)	3.90					
Operating Control Conditions	TX Level	(dB)	0.00	0.00		0.00		0.00
	Focus	(mm)	15.00	5.00		5.00		5.00
a	This Index is not relevant to this operating mode.							
b	This transducer is not intended for transcranial or neonatal cephalic uses.							
c	This formulation for TIS is less than that for an alternate formulation in this mode.							
#	No data is provided for this operation condition since the maximum index value is not reported for the reason listed.							

Acoustic Output Reporting Table — Track 3, FDA 510(k) and IEC 60601-2-37

(Per transducer/mode that exceeds MI or TI value of 1.0)

Transducer model: **9MC3**Operating mode: **M-mode**

Index Label			MI	TIS		TIB		TIC
				At surface	Below surface	At surface	Below surface	
Maximum Index Value			1.50	0.85		1.19		0.89
Index Component Value				0.85	0.46	0.75	1.19	
Acoustic Parameters	IEC	Units						
	$p_{ra,\alpha}$ at z_{MI}	(MPa)	3.53					
	P	(mW)		#		35.08		35.08
	P_{1x1}	(mW)		32.45		35.08		
	z_s	(cm)			1.60			
	z_b	(cm)					2.10	
	z_{MI}	(cm)	1.70					
	$z_{pii,\alpha}$	(cm)	1.70					
	f_{awf}	(MHz)	5.56	5.53		4.49		4.49
Other Information	p_{rr}	(Hz)	720.00					
	s_{rr}	(Hz)	#					
	n_{pps}	-	#					
	$I_{pa,\alpha}$ at $z_{pii,\alpha}$	(W/cm ²)	549.43					
	$I_{spta,\alpha}$ at $z_{pii,\alpha}$ or $z_{sii,\alpha}$	(mW/cm ²)	176.37					
	I_{spta} at z_{pii} or z_{sii}	(mW/cm ²)	339.36					
	p_r at z_{pii}	(MPa)	4.64					
Operating Control Conditions	TX Level	(dB)	0.00	0.00		0.00		0.00
	Focus	(mm)	20.00	80.00		40.00		40.00
a This Index is not relevant to this operating mode.								
b This transducer is not intended for transcranial or neonatal cephalic uses.								
c This formulation for TIS is less than that for an alternate formulation in this mode.								
# No data is provided for this operation condition since the maximum index value is not reported for the reason listed.								

Acoustic Output Reporting Table — Track 3, FDA 510(k) and IEC 60601-2-37

(Per transducer/mode that exceeds MI or TI value of 1.0)

Transducer model: **9MC3**Operating mode: **Pulsed Wave Doppler**

Index Label			MI	TIS		TIB		TIC
				At surface	Below surface	At surface	Below surface	
Maximum Index Value			1.48	1.00		2.07		1.50
Index Component Value				1.00	0.55	0.55	2.07	
Acoustic Parameters	IEC	Units						
	$p_{ra,\alpha}$ at z_{MI}	(MPa)	2.99					
	P	(mW)		#		28.08		24.23
	$P_{1\times1}$	(mW)		41.35		28.08		
	z_s	(cm)			1.80			
	z_b	(cm)					0.90	
	z_{MI}	(cm)	1.40					
	$z_{pii,\alpha}$	(cm)	1.40					
	f_{awf}	(MHz)	4.08	5.10		4.09		5.29
Other Information	p_{rr}	(Hz)	976.56					
	s_{rr}	(Hz)	#					
	η_{pps}	-	#					
	$I_{pa,\alpha}$ at $z_{pii,\alpha}$	(W/cm ²)	544.67					
	$I_{spta,\alpha}$ at $z_{pii,\alpha}$ or $z_{sii,\alpha}$	(mW/cm ²)	496.24					
	I_{spta} at z_{pii} or z_{sii}	(mW/cm ²)	768.16					
	p_r at z_{pii}	(MPa)	3.60					
	Operating Control Conditions	TX Level	(dB)	0.00	0.00		0.00	
Focus		(mm)	20.00	60.00		10.00		5.00
a This Index is not relevant to this operating mode.								
b This transducer is not intended for transcranial or neonatal cephalic uses.								
c This formulation for TIS is less than that for an alternate formulation in this mode.								
# No data is provided for this operation condition since the maximum index value is not reported for the reason listed.								

Acoustic Output Reporting Table — Track 3, FDA 510(k) and IEC 60601-2-37

(Per transducer/mode that exceeds MI or TI value of 1.0)

Transducer model: **9MC3**Operating mode: **Color Flow Doppler**

Index Label			MI	TIS		TIB		TIC
				At surface	Below surface	At surface	Below surface	
Maximum Index Value			1.59	0.43		0.43		0.43
Index Component Value				0.43	0.43	0.43	0.43	
Acoustic Parameters	IEC	Units						
	$p_{ra,\alpha}$ at z_{MI}	(MPa)	3.62					
	P	(mW)		#		18.83		22.68
	$P_{1\times1}$	(mW)		17.61		17.61		
	z_s	(cm)			#			
	z_b	(cm)					#	
	z_{MI}	(cm)	1.00					
	$z_{pII,\alpha}$	(cm)	1.00					
	f_{awf}	(MHz)	5.22	5.10		5.10		5.24
Other Information	p_{rr}	(Hz)	#					
	s_{rr}	(Hz)	65.75					
	n_{pps}	-	10.00					
	$I_{pa,\alpha}$ at $z_{pII,\alpha}$	(W/cm ²)	636.22					
	$I_{spta,\alpha}$ at $z_{pII,\alpha}$ or $z_{sII,\alpha}$	(mW/cm ²)	23.55					
	I_{spta} at z_{pII} or z_{sII}	(mW/cm ²)	34.99					
	p_r at z_{pII}	(MPa)	4.18					
Operating Control Conditions	TX Level	(dB)	0.00	0.00		0.00		0.00
	Focus	(mm)	15.00	30.00		30.00		10.00
a	This Index is not relevant to this operating mode.							
b	This transducer is not intended for transcranial or neonatal cephalic uses.							
c	This formulation for TIS is less than that for an alternate formulation in this mode.							
#	No data is provided for this operation condition since the maximum index value is not reported for the reason listed.							

Acoustic Output Reporting Table — Track 3, FDA 510(k) and IEC 60601-2-37

(Per transducer/mode that exceeds MI or TI value of 1.0)

Transducer model: **10MC3**Operating mode: **2D-mode**

Index Label			MI	TIS		TIB		TIC
				At surface	Below surface	At surface	Below surface	
Maximum Index Value			1.11	0.73		0.73		0.73
Index Component Value				0.73	0.73	0.73	0.73	
Acoustic Parameters	IEC	Units						
	$p_{ra,\alpha}$ at z_{MI}	(MPa)	2.29					
	P	(mW)		#		24.31		41.53
	$P_{1\times1}$	(mW)		22.34		22.34		
	z_s	(cm)			#			
	z_b	(cm)					#	
	z_{MI}	(cm)	1.90					
	$z_{pii,\alpha}$	(cm)	1.90					
	f_{awf}	(MHz)	4.37	6.86		6.86		6.07
Other Information	prr	(Hz)	#					
	srr	(Hz)	49.96					
	n_{pps}	-	2.00					
	$I_{pa,\alpha}$ at $z_{pii,\alpha}$	(W/cm ²)	224.13					
	$I_{spta,\alpha}$ at $z_{pii,\alpha}$ or $z_{sii,\alpha}$	(mW/cm ²)	8.09					
	I_{spta} at z_{pii} or z_{sii}	(mW/cm ²)	13.10					
	p_r at z_{pii}	(MPa)	2.96					
Operating Control Conditions	TX Level	(dB)	0.00	0.00		0.00		0.00
	Focus	(mm)	20.00	15.00		15.00		10.00
a	This Index is not relevant to this operating mode.							
b	This transducer is not intended for transcranial or neonatal cephalic uses.							
c	This formulation for TIS is less than that for an alternate formulation in this mode.							
#	No data is provided for this operation condition since the maximum index value is not reported for the reason listed.							

Acoustic Output Reporting Table — Track 3, FDA 510(k) and IEC 60601-2-37

(Per transducer/mode that exceeds MI or TI value of 1.0)

Transducer model: **10MC3**Operating mode: **M-mode**

Index Label			MI	TIS		TIB		TIC
				At surface	Below surface	At surface	Below surface	
Maximum Index Value			1.43	0.57		0.73		0.74
Index Component Value				0.57	0.36	0.32	0.73	
Acoustic Parameters	IEC	Units						
	$p_{ra,\alpha}$ at z_{MI}	(MPa)	3.25					
	P	(mW)		#		13.32		21.55
	P_{1x1}	(mW)		19.80		13.32		
	z_s	(cm)			1.10			
	z_b	(cm)					0.90	
	z_{MI}	(cm)	0.80					
	$z_{pii,\alpha}$	(cm)	0.80					
	f_{awf}	(MHz)	5.02	6.05		5.02		4.82
Other Information	p_{rr}	(Hz)	720.00					
	s_{rr}	(Hz)	#					
	n_{pps}	-	#					
	$I_{pa,\alpha}$ at $z_{pii,\alpha}$	(W/cm ²)	425.58					
	$I_{spta,\alpha}$ at $z_{pii,\alpha}$ or $z_{sii,\alpha}$	(mW/cm ²)	136.78					
	I_{spta} at z_{pii} or z_{sii}	(mW/cm ²)	161.74					
	p_r at z_{pii}	(MPa)	3.35					
Operating Control Conditions	TX Level	(dB)	0.00	0.00		0.00		0.00
	Focus	(mm)	15.00	80.00		15.00		80.00
a This Index is not relevant to this operating mode.								
b This transducer is not intended for transcranial or neonatal cephalic uses.								
c This formulation for TIS is less than that for an alternate formulation in this mode.								
# No data is provided for this operation condition since the maximum index value is not reported for the reason listed.								

Acoustic Output Reporting Table — Track 3, FDA 510(k) and IEC 60601-2-37

(Per transducer/mode that exceeds MI or TI value of 1.0)

Transducer model: **10MC3**Operating mode: **Pulsed Wave Doppler**

Index Label			MI	TIS		TIB		TIC
				At surface	Below surface	At surface	Below surface	
Maximum Index Value			1.50	0.74		1.23		1.32
Index Component Value				0.74	0.41	0.37	1.23	
Acoustic Parameters	IEC	Units						
	$p_{ra,\alpha}$ at z_{MI}	(MPa)	3.03					
	P	(mW)		#		19.08		16.35
	P_{1x1}	(mW)		26.03		19.08		
	z_s	(cm)			1.40			
	z_b	(cm)					1.60	
	z_{MI}	(cm)	1.70					
	$z_{pII,\alpha}$	(cm)	1.70					
	f_{awf}	(MHz)	4.08	5.93		4.08		5.29
	prf	(Hz)	610.35					
Other Information	srf	(Hz)	#					
	n_{pps}	-	#					
	$I_{pa,\alpha}$ at $z_{pII,\alpha}$	(W/cm ²)	547.59					
	$I_{spta,\alpha}$ at $z_{pII,\alpha}$ or $z_{sII,\alpha}$	(mW/cm ²)	243.70					
	I_{spta} at z_{pII} or z_{sII}	(mW/cm ²)	393.50					
	p_r at z_{pII}	(MPa)	3.85					
	TX Level	(dB)	0.00	0.00		0.00		0.00
	Focus	(mm)	20.00	80.00		20.00		5.00
Operating Control Conditions								
a This Index is not relevant to this operating mode. b This transducer is not intended for transcranial or neonatal cephalic uses. c This formulation for TIS is less than that for an alternate formulation in this mode. # No data is provided for this operation condition since the maximum index value is not reported for the reason listed.								

Acoustic Output Reporting Table — Track 3, FDA 510(k) and IEC 60601-2-37

(Per transducer/mode that exceeds MI or TI value of 1.0)

Transducer model: **10MC3**Operating mode: **Color Flow Doppler**

Index Label			MI	TIS		TIB		TIC
				At surface	Below surface	At surface	Below surface	
Maximum Index Value			1.34	0.34		0.34		0.66
Index Component Value				0.34	0.34	0.34	0.34	
Acoustic Parameters	IEC	Units						
	$p_{ra,\alpha}$ at z_{MI}	(MPa)	3.07					
	P	(mW)		#		13.28		37.71
	P_{1x1}	(mW)		12.05		12.05		
	z_s	(cm)			#			
	z_b	(cm)					#	
	z_{MI}	(cm)	0.60					
	$z_{pII,\alpha}$	(cm)	0.70					
	f_{awf}	(MHz)	5.26	6.00		6.00		5.26
	prf	(Hz)	#					
Other Information	srf	(Hz)	40.01					
	n_{pps}	-	10.00					
	$I_{pa,\alpha}$ at $z_{pII,\alpha}$	(W/cm ²)	408.24					
	$I_{spta,\alpha}$ at $z_{pII,\alpha}$ or $z_{sII,\alpha}$	(mW/cm ²)	11.63					
	I_{spta} at z_{pII} or z_{sII}	(mW/cm ²)	15.55					
	p_r at z_{pII}	(MPa)	3.28					
	TX Level	(dB)	0.00	0.00		0.00		0.00
	Focus	(mm)	10.00	10.00		10.00		10.00
Operating Control Conditions								
a This Index is not relevant to this operating mode. b This transducer is not intended for transcranial or neonatal cephalic uses. c This formulation for TIS is less than that for an alternate formulation in this mode. # No data is provided for this operation condition since the maximum index value is not reported for the reason listed.								

Acoustic Output Reporting Table — Track 3, FDA 510(k) and IEC 60601-2-37

(Per transducer/mode that exceeds MI or TI value of 1.0)

Transducer model: **10BCC3**Operating mode: **M-mode**

Index Label			MI	TIS		TIB		TIC
				At surface	Below surface	At surface	Below surface	
Maximum Index Value			1.09	0.30		0.78		0.39
Index Component Value				0.30	0.19	0.21	0.78	
Acoustic Parameters	IEC	Units						
	$p_{ra,\alpha}$ at z_{MI}	(MPa)	2.06					
	P	(mW)		#		11.50		12.28
	P_{1x1}	(mW)		14.65		#		
	z_s	(cm)			1.50			
	z_b	(cm)					1.80	
	z_{MI}	(cm)	1.90					
	$z_{pii,\alpha}$	(cm)	1.90					
	f_{awf}	(MHz)	3.69	4.31		3.92		4.28
Other Information	p_{rr}	(Hz)	1440.09					
	s_{rr}	(Hz)	#					
	η_{pps}	-	#					
	$I_{pa,\alpha}$ at $z_{pii,\alpha}$	(W/cm ²)	234.49					
	$I_{spta,\alpha}$ at $z_{pii,\alpha}$ or $z_{sii,\alpha}$	(mW/cm ²)	179.12					
	I_{spta} at z_{pii} or z_{sii}	(mW/cm ²)	293.53					
	p_r at z_{pii}	(MPa)	2.64					
Operating Control Conditions	TX Level	(dB)	0.00	0.00		0.00		0.00
	Focus	(mm)	20.00	80.00		20.00		20.00
a This Index is not relevant to this operating mode.								
b This transducer is not intended for transcranial or neonatal cephalic uses.								
c This formulation for TIS is less than that for an alternate formulation in this mode.								
# No data is provided for this operation condition since the maximum index value is not reported for the reason listed.								

Acoustic Output Reporting Table — Track 3, FDA 510(k) and IEC 60601-2-37

(Per transducer/mode that exceeds MI or TI value of 1.0)

Transducer model: **10BCC3**Operating mode: **Pulsed Wave Doppler**

Index Label			MI	TIS		TIB		TIC
				At surface	Below surface	At surface	Below surface	
Maximum Index Value			1.17	0.39		1.20		0.93
Index Component Value				0.39	0.27	0.35	1.20	
Acoustic Parameters	IEC	Units						
	$p_{ra,\alpha}$ at z_{MI}	(MPa)	2.46					
	P	(mW)		#		16.63		12.11
	P_{1x1}	(mW)		18.14		#		
	z_s	(cm)			1.20			
	z_b	(cm)					1.70	
	z_{MI}	(cm)	2.10					
	$z_{pii,\alpha}$	(cm)	2.10					
	f_{awf}	(MHz)	4.46	4.46		4.47		4.46
Other Information	p_{rr}	(Hz)	390.63					
	s_{rr}	(Hz)	#					
	n_{pps}	-	#					
	$I_{pa,\alpha}$ at $z_{pii,\alpha}$	(W/cm ²)	419.58					
	$I_{spta,\alpha}$ at $z_{pii,\alpha}$ or $z_{sii,\alpha}$	(mW/cm ²)	141.67					
	I_{spta} at z_{pii} or z_{sii}	(mW/cm ²)	262.35					
	p_r at z_{pii}	(MPa)	3.27					
Operating Control Conditions	TX Level	(dB)	0.00	0.00		0.00		0.00
	Focus	(mm)	30.00	30.00		20.00		5.00
a	This Index is not relevant to this operating mode.							
b	This transducer is not intended for transcranial or neonatal cephalic uses.							
c	This formulation for TIS is less than that for an alternate formulation in this mode.							
#	No data is provided for this operation condition since the maximum index value is not reported for the reason listed.							

Acoustic Output Reporting Table — Track 3, FDA 510(k) and IEC 60601-2-37

(Per transducer/mode that exceeds MI or TI value of 1.0)

Transducer model: **10BCC3**Operating mode: **Color Flow Doppler**

Index Label			MI	TIS		TIB		TIC
				At surface	Below surface	At surface	Below surface	
Maximum Index Value			1.07	0.33		0.33		0.45
Index Component Value				0.29	0.29	0.29	0.29	
Acoustic Parameters	IEC	Units						
	$p_{ra,\alpha}$ at z_{MI}	(MPa)	2.37					
	P	(mW)		#		13.66		22.18
	$P_{1\times1}$	(mW)		12.63		12.63		
	z_s	(cm)			#			
	z_b	(cm)					#	
	z_{MI}	(cm)	1.90					
	$z_{pii,\alpha}$	(cm)	1.90					
	f_{awf}	(MHz)	4.92	4.89		4.89		4.43
Other Information	prr	(Hz)	#					
	srr	(Hz)	75.09					
	n_{pps}	-	14.00					
	$I_{pa,\alpha}$ at $z_{pii,\alpha}$	(W/cm ²)	357.54					
	$I_{spta,\alpha}$ at $z_{pii,\alpha}$ or $z_{sii,\alpha}$	(mW/cm ²)	11.43					
	I_{spta} at z_{pii} or z_{sii}	(mW/cm ²)	21.80					
	p_r at z_{pii}	(MPa)	3.16					
Operating Control Conditions	TX Level	(dB)	0.00	0.00		0.00		0.00
	Focus	(mm)	20.00	80.00		80.00		80.00
a	This Index is not relevant to this operating mode.							
b	This transducer is not intended for transcranial or neonatal cephalic uses.							
c	This formulation for TIS is less than that for an alternate formulation in this mode.							
#	No data is provided for this operation condition since the maximum index value is not reported for the reason listed.							

Acoustic Output Reporting Table — Track 3, FDA 510(k) and IEC 60601-2-37

(Per transducer/mode that exceeds MI or TI value of 1.0)

Transducer model: **9VC2**Operating mode: **2D-mode**

Index Label			MI	TIS		TIB		TIC
				At surface	Below surface	At surface	Below surface	
Maximum Index Value			1.40	0.89		0.89		1.33
Index Component Value				0.89	0.89	0.89	0.89	
Acoustic Parameters	IEC	Units						
	$p_{ra,\alpha}$ at z_{MI}	(MPa)	2.18					
	P	(mW)		#		48.48		139.15
	P_{1x1}	(mW)		40.56		40.56		
	z_s	(cm)			#			
	z_b	(cm)					#	
	z_{MI}	(cm)	3.20					
	$z_{pii,\alpha}$	(cm)	3.20					
	f_{awf}	(MHz)	2.42	4.60		4.60		2.94
Other Information	p_{rr}	(Hz)	#					
	s_{rr}	(Hz)	35.27					
	n_{pps}	-	1.00					
	$I_{pa,\alpha}$ at $z_{pii,\alpha}$	(W/cm ²)	322.37					
	$I_{spta,\alpha}$ at $z_{pii,\alpha}$ or $z_{sii,\alpha}$	(mW/cm ²)	12.30					
	I_{spta} at z_{pii} or z_{sii}	(mW/cm ²)	20.97					
	p_r at z_{pii}	(MPa)	2.86					
Operating Control Conditions	TX Level	(dB)	0.00	0.00		0.00		0.00
	Focus	(mm)	40.00	20.00		20.00		200.00
a This Index is not relevant to this operating mode.								
b This transducer is not intended for transcranial or neonatal cephalic uses.								
c This formulation for TIS is less than that for an alternate formulation in this mode.								
# No data is provided for this operation condition since the maximum index value is not reported for the reason listed.								

Acoustic Output Reporting Table — Track 3, FDA 510(k) and IEC 60601-2-37

(Per transducer/mode that exceeds MI or TI value of 1.0)

Transducer model: **9VC2**Operating mode: **M-mode**

Index Label			MI	TIS		TIB		TIC
				At surface	Below surface	At surface	Below surface	
Maximum Index Value			1.39	0.69		1.52		1.13
Index Component Value				0.33	0.69	0.43	1.52	
Acoustic Parameters	IEC	Units						
	$p_{ra,\alpha}$ at z_{MI}	(MPa)	2.13					
	P	(mW)		#		62.44		109.15
	P_{1x1}	(mW)		25.13		33.19		
	z_s	(cm)			3.50			
	z_b	(cm)					4.20	
	z_{MI}	(cm)	4.20					
	$z_{pii,\alpha}$	(cm)	4.20					
	f_{awf}	(MHz)	2.38	2.80		2.74		3.39
Other Information	pr_r	(Hz)	720.00					
	srr	(Hz)	#					
	n_{pps}	-	#					
	$I_{pa,\alpha}$ at $z_{pii,\alpha}$	(W/cm ²)	347.15					
	$I_{spta,\alpha}$ at $z_{pii,\alpha}$ or $z_{sii,\alpha}$	(mW/cm ²)	174.10					
	I_{spta} at z_{pii} or z_{sii}	(mW/cm ²)	372.13					
	p_r at z_{pii}	(MPa)	3.07					
	Operating Control Conditions	TX Level	(dB)	0.00	0.00		0.00	
Focus		(mm)	60.00	160.00		60.00		200.00
a This Index is not relevant to this operating mode.								
b This transducer is not intended for transcranial or neonatal cephalic uses.								
c This formulation for TIS is less than that for an alternate formulation in this mode.								
# No data is provided for this operation condition since the maximum index value is not reported for the reason listed.								

Acoustic Output Reporting Table — Track 3, FDA 510(k) and IEC 60601-2-37

(Per transducer/mode that exceeds MI or TI value of 1.0)

Transducer model: **9VC2**Operating mode: **Pulsed Wave Doppler**

Index Label			MI	TIS		TIB		TIC
				At surface	Below surface	At surface	Below surface	
Maximum Index Value			1.40	0.76		1.86		1.23
Index Component Value				0.76	0.54	0.61	1.86	
Acoustic Parameters	IEC	Units						
	$p_{ra,\alpha}$ at z_{MI}	(MPa)	2.50					
	P	(mW)		#		50.64		37.66
	P_{1x1}	(mW)		49.76		41.77		
	z_s	(cm)			1.60			
	z_b	(cm)					3.20	
	z_{MI}	(cm)	3.60					
	$z_{pii,\alpha}$	(cm)	3.60					
	f_{awf}	(MHz)	3.23	3.22		3.04		3.27
Other Information	p_{rr}	(Hz)	610.35					
	s_{rr}	(Hz)	#					
	n_{pps}	-	#					
	$I_{pa,\alpha}$ at $z_{pii,\alpha}$	(W/cm ²)	349.71					
	$I_{spta,\alpha}$ at $z_{pii,\alpha}$ or $z_{sii,\alpha}$	(mW/cm ²)	126.43					
	I_{spta} at z_{pii} or z_{sii}	(mW/cm ²)	285.93					
	p_r at z_{pii}	(MPa)	3.51					
Operating Control Conditions	TX Level	(dB)	0.00	0.00		0.00		0.00
	Focus	(mm)	60.00	40.00		55.00		20.00
a	This Index is not relevant to this operating mode.							
b	This transducer is not intended for transcranial or neonatal cephalic uses.							
c	This formulation for TIS is less than that for an alternate formulation in this mode.							
#	No data is provided for this operation condition since the maximum index value is not reported for the reason listed.							

Acoustic Output Reporting Table — Track 3, FDA 510(k) and IEC 60601-2-37

(Per transducer/mode that exceeds MI or TI value of 1.0)

Transducer model: **9VC2**Operating mode: **Color Flow Doppler**

Index Label			MI	TIS		TIB		TIC
				At surface	Below surface	At surface	Below surface	
Maximum Index Value			1.40	0.73		0.73		1.33
Index Component Value				0.31	0.31	0.31	0.31	
Acoustic Parameters	IEC	Units						
	$p_{ra,\alpha}$ at z_{MI}	(MPa)	2.50					
	P	(mW)		#		23.90		67.14
	P_{1x1}	(mW)		19.58		19.58		
	z_s	(cm)			#			
	z_b	(cm)					#	
	z_{MI}	(cm)	0.50					
	$z_{pii,\alpha}$	(cm)	1.10					
	f_{awf}	(MHz)	3.27	3.27		3.27		3.03
Other Information	pr_r	(Hz)	#					
	sr_r	(Hz)	68.10					
	n_{pps}	-	10.00					
	$I_{pa,\alpha}$ at $z_{pii,\alpha}$	(W/cm ²)	220.53					
	$I_{spta,\alpha}$ at $z_{pii,\alpha}$ or $z_{sii,\alpha}$	(mW/cm ²)	9.52					
	I_{spta} at z_{pii} or z_{sii}	(mW/cm ²)	21.95					
	p_r at z_{pii}	(MPa)	2.39					
Operating Control Conditions	TX Level	(dB)	0.00	0.00		0.00		0.00
	Focus	(mm)	20.00	20.00		20.00		40.00
a This Index is not relevant to this operating mode.								
b This transducer is not intended for transcranial or neonatal cephalic uses.								
c This formulation for TIS is less than that for an alternate formulation in this mode.								
# No data is provided for this operation condition since the maximum index value is not reported for the reason listed.								

Acoustic Output Reporting Table — Track 3, FDA 510(k) and IEC 60601-2-37

(Per transducer/mode that exceeds MI or TI value of 1.0)

Transducer model: **8VC3**Operating mode: **2D-mode**

Index Label			MI	TIS		TIB		TIC
				At surface	Below surface	At surface	Below surface	
Maximum Index Value			1.49	1.76		1.76		3.42
Index Component Value				1.76	1.76	1.76	1.76	
Acoustic Parameters	IEC	Units						
	$p_{ra,\alpha}$ at z_{MI}	(MPa)	2.31					
	P	(mW)		#		142.34		369.47
	P_{1x1}	(mW)		102.12		102.12		
	z_s	(cm)			#			
	z_b	(cm)					#	
	z_{MI}	(cm)	3.80					
	$z_{pii,\alpha}$	(cm)	3.80					
	f_{awf}	(MHz)	2.41	3.62		3.62		2.82
Other Information	p_{rr}	(Hz)	#					
	s_{rr}	(Hz)	21.73					
	η_{pps}	-	2.00					
	$I_{pa,\alpha}$ at $z_{pii,\alpha}$	(W/cm ²)	269.53					
	$I_{spta,\alpha}$ at $z_{pii,\alpha}$ or $z_{sii,\alpha}$	(mW/cm ²)	9.21					
	I_{spta} at z_{pii} or z_{sii}	(mW/cm ²)	17.91					
	p_r at z_{pii}	(MPa)	3.06					
Operating Control Conditions	TX Level	(dB)	0.00	0.00		0.00		0.00
	Focus	(mm)	40.00	20.00		20.00		80.00
a	This Index is not relevant to this operating mode.							
b	This transducer is not intended for transcranial or neonatal cephalic uses.							
c	This formulation for TIS is less than that for an alternate formulation in this mode.							
#	No data is provided for this operation condition since the maximum index value is not reported for the reason listed.							

Acoustic Output Reporting Table — Track 3, FDA 510(k) and IEC 60601-2-37

(Per transducer/mode that exceeds MI or TI value of 1.0)

Transducer model: **8VC3**Operating mode: **M-mode**

Index Label			MI	TIS		TIB		TIC
				At surface	Below surface	At surface	Below surface	
Maximum Index Value			1.49	1.19		2.24		2.49
Index Component Value				0.91	1.19	0.91	2.24	
Acoustic Parameters	IEC	Units						
	$p_{ra,\alpha}$ at z_{MI}	(MPa)	2.31					
	P	(mW)		#		117.28		213.00
	$P_{1\times1}$	(mW)		67.68		67.68		
	z_s	(cm)			3.60			
	z_b	(cm)					3.80	
	z_{MI}	(cm)	3.80					
	$z_{pii,\alpha}$	(cm)	3.80					
	f_{awf}	(MHz)	2.41	2.84		2.84		2.82
Other Information	prr	(Hz)	720.00					
	srr	(Hz)	#					
	n_{pps}	-	#					
	$I_{pa,\alpha}$ at $z_{pii,\alpha}$	(W/cm ²)	269.53					
	$I_{spta,\alpha}$ at $z_{pii,\alpha}$ or $z_{sii,\alpha}$	(mW/cm ²)	168.89					
	I_{spta} at z_{pii} or z_{sii}	(mW/cm ²)	314.57					
	p_r at z_{pii}	(MPa)	3.06					
Operating Control Conditions	TX Level	(dB)	0.00	0.00		0.00		0.00
	Focus	(mm)	40.00	60.00		40.00		80.00
a	This Index is not relevant to this operating mode.							
b	This transducer is not intended for transcranial or neonatal cephalic uses.							
c	This formulation for TIS is less than that for an alternate formulation in this mode.							
#	No data is provided for this operation condition since the maximum index value is not reported for the reason listed.							

Acoustic Output Reporting Table — Track 3, FDA 510(k) and IEC 60601-2-37

(Per transducer/mode that exceeds MI or TI value of 1.0)

Transducer model: **8VC3**Operating mode: **Pulsed Wave Doppler**

Index Label			MI	TIS		TIB		TIC
				At surface	Below surface	At surface	Below surface	
Maximum Index Value			1.50	0.94		2.00		1.54
Index Component Value				0.94	0.75	0.67	2.00	
Acoustic Parameters	IEC	Units						
	$p_{ra,\alpha}$ at z_{MI}	(MPa)	2.45					
	P	(mW)		#		53.04		53.04
	$P_{1\times1}$	(mW)		61.44		53.04		
	z_s	(cm)			1.90			
	z_b	(cm)					1.30	
	z_{MI}	(cm)	4.20					
	$z_{pii,\alpha}$	(cm)	4.20					
	f_{awf}	(MHz)	2.67	3.22		2.67		2.67
Other Information	p_{rr}	(Hz)	390.63					
	s_{rr}	(Hz)	#					
	n_{pps}	-	#					
	$I_{pa,\alpha}$ at $z_{pii,\alpha}$	(W/cm ²)	430.45					
	$I_{spta,\alpha}$ at $z_{pii,\alpha}$ or $z_{sii,\alpha}$	(mW/cm ²)	248.65					
	I_{spta} at z_{pii} or z_{sii}	(mW/cm ²)	510.25					
	p_r at z_{pii}	(MPa)	3.48					
Operating Control Conditions	TX Level	(dB)	0.00	0.00		0.00		0.00
	Focus	(mm)	60.00	40.00		20.00		20.00
a This Index is not relevant to this operating mode.								
b This transducer is not intended for transcranial or neonatal cephalic uses.								
c This formulation for TIS is less than that for an alternate formulation in this mode.								
# No data is provided for this operation condition since the maximum index value is not reported for the reason listed.								

Acoustic Output Reporting Table — Track 3, FDA 510(k) and IEC 60601-2-37

(Per transducer/mode that exceeds MI or TI value of 1.0)

Transducer model: **8VC3**Operating mode: **Color Flow Doppler**

Index Label			MI	TIS		TIB		TIC
				At surface	Below surface	At surface	Below surface	
Maximum Index Value			1.50	0.39		0.39		0.97
Index Component Value				0.39	0.39	0.39	0.39	
Acoustic Parameters	IEC	Units						
	$p_{ra,\alpha}$ at z_{MI}	(MPa)	2.44					
	P	(mW)		#		40.14		105.42
	$P_{1\times1}$	(mW)		25.31		25.31		
	z_s	(cm)			#			
	z_b	(cm)					#	
	z_{MI}	(cm)	3.50					
	$z_{pii,\alpha}$	(cm)	3.50					
	f_{awf}	(MHz)	2.68	3.22		3.22		2.67
Other Information	prr	(Hz)	#					
	srr	(Hz)	70.19					
	n_{pps}	-	10.00					
	$I_{pa,\alpha}$ at $z_{pii,\alpha}$	(W/cm ²)	401.91					
	$I_{spta,\alpha}$ at $z_{pii,\alpha}$ or $z_{sii,\alpha}$	(mW/cm ²)	25.93					
	I_{spta} at z_{pii} or z_{sii}	(mW/cm ²)	53.30					
	p_r at z_{pii}	(MPa)	3.20					
Operating Control Conditions	TX Level	(dB)	0.00	0.00		0.00		0.00
	Focus	(mm)	40.00	40.00		40.00		200.00
a	This Index is not relevant to this operating mode.							
b	This transducer is not intended for transcranial or neonatal cephalic uses.							
c	This formulation for TIS is less than that for an alternate formulation in this mode.							
#	No data is provided for this operation condition since the maximum index value is not reported for the reason listed.							

Acoustic Output Reporting Table — Track 3, FDA 510(k) and IEC 60601-2-37

(Per transducer/mode that exceeds MI or TI value of 1.0)

Transducer model: **9VE3**Operating mode: **2D-mode**

Index Label			MI	TIS		TIB		TIC
				At surface	Below surface	At surface	Below surface	
Maximum Index Value			1.41	0.64		0.64		0.72
Index Component Value				0.64	0.64	0.64	0.64	
Acoustic Parameters	IEC	Units						
	$p_{ra,\alpha}$ at z_{MI}	(MPa)	3.24					
	P	(mW)		#		25.82		28.56
	P_{1x1}	(mW)		24.75		24.75		
	z_s	(cm)			#			
	z_b	(cm)					#	
	z_{MI}	(cm)	0.60					
	$z_{pii,\alpha}$	(cm)	0.70					
	f_{awf}	(MHz)	5.47	5.47		5.47		4.61
Other Information	p_{rr}	(Hz)	#					
	s_{rr}	(Hz)	17.34					
	n_{pps}	-	1.00					
	$I_{pa,\alpha}$ at $z_{pii,\alpha}$	(W/cm ²)	424.86					
	$I_{spta,\alpha}$ at $z_{pii,\alpha}$ or $z_{sii,\alpha}$	(mW/cm ²)	30.45					
	I_{spta} at z_{pii} or z_{sii}	(mW/cm ²)	41.18					
	p_r at z_{pii}	(MPa)	3.47					
Operating Control Conditions	TX Level	(dB)	0.00	0.00		0.00		0.00
	Focus	(mm)	10.00	10.00		10.00		20.00
a This Index is not relevant to this operating mode.								
b This transducer is not intended for transcranial or neonatal cephalic uses.								
c This formulation for TIS is less than that for an alternate formulation in this mode.								
# No data is provided for this operation condition since the maximum index value is not reported for the reason listed.								

Acoustic Output Reporting Table — Track 3, FDA 510(k) and IEC 60601-2-37

(Per transducer/mode that exceeds MI or TI value of 1.0)

Transducer model: **9VE3**Operating mode: **M-mode**

Index Label			MI	TIS		TIB		TIC
				At surface	Below surface	At surface	Below surface	
Maximum Index Value			1.50	0.65		0.92		0.81
Index Component Value				0.65	0.41	0.63	0.92	
Acoustic Parameters	IEC	Units						
	$p_{ra,\alpha}$ at z_{MI}	(MPa)	3.35					
	P	(mW)		#		28.81		28.81
	$P_{1\times1}$	(mW)		29.64		28.81		
	z_s	(cm)			1.50			
	z_b	(cm)					1.90	
	z_{MI}	(cm)	1.80					
	$z_{pii,\alpha}$	(cm)	1.80					
	f_{awf}	(MHz)	5.06	4.64		4.56		4.56
Other Information	prr	(Hz)	720.00					
	srr	(Hz)	#					
	n_{pps}	-	#					
	$I_{pa,\alpha}$ at $z_{pii,\alpha}$	(W/cm ²)	499.33					
	$I_{spta,\alpha}$ at $z_{pii,\alpha}$ or $z_{sii,\alpha}$	(mW/cm ²)	107.07					
	I_{spta} at z_{pii} or z_{sii}	(mW/cm ²)	196.56					
	p_r at z_{pii}	(MPa)	4.27					
	Operating Control Conditions	TX Level	(dB)	0.00	0.00		0.00	
Focus		(mm)	30.00	60.00		40.00		40.00
a	This Index is not relevant to this operating mode.							
b	This transducer is not intended for transcranial or neonatal cephalic uses.							
c	This formulation for TIS is less than that for an alternate formulation in this mode.							
#	No data is provided for this operation condition since the maximum index value is not reported for the reason listed.							

Acoustic Output Reporting Table — Track 3, FDA 510(k) and IEC 60601-2-37

(Per transducer/mode that exceeds MI or TI value of 1.0)

Transducer model: **9VE3**Operating mode: **Pulsed Wave Doppler**

Index Label			MI	TIS		TIB		TIC
				At surface	Below surface	At surface	Below surface	
Maximum Index Value			1.50	0.39		1.01		0.73
Index Component Value				0.39	0.25	0.32	1.01	
Acoustic Parameters	IEC	Units						
	$p_{ra,\alpha}$ at z_{MI}	(MPa)	2.99					
	P	(mW)		#		13.20		8.74
	P_{1x1}	(mW)		16.21		13.20		
	z_s	(cm)			1.30			
	z_b	(cm)					0.90	
	z_{MI}	(cm)	1.30					
	$z_{pii,\alpha}$	(cm)	1.30					
	f_{awf}	(MHz)	3.98	5.08		5.08		5.15
Other Information	p_{rr}	(Hz)	151.99					
	s_{rr}	(Hz)	#					
	n_{pps}	-	#					
	$I_{pa,\alpha}$ at $z_{pii,\alpha}$	(W/cm ²)	511.97					
	$I_{spta,\alpha}$ at $z_{pii,\alpha}$ or $z_{sii,\alpha}$	(mW/cm ²)	54.36					
	I_{spta} at z_{pii} or z_{sii}	(mW/cm ²)	76.69					
	p_r at z_{pii}	(MPa)	3.39					
Operating Control Conditions	TX Level	(dB)	0.00	0.00		0.00		0.00
	Focus	(mm)	15.00	80.00		15.00		5.00
a This Index is not relevant to this operating mode.								
b This transducer is not intended for transcranial or neonatal cephalic uses.								
c This formulation for TIS is less than that for an alternate formulation in this mode.								
# No data is provided for this operation condition since the maximum index value is not reported for the reason listed.								

Acoustic Output Reporting Table — Track 3, FDA 510(k) and IEC 60601-2-37

(Per transducer/mode that exceeds MI or TI value of 1.0)

Transducer model: **9VE3**Operating mode: **Color Flow Doppler**

Index Label			MI	TIS		TIB		TIC
				At surface	Below surface	At surface	Below surface	
Maximum Index Value			1.50	0.19		0.19		0.21
Index Component Value				0.19	0.19	0.19	0.19	
Acoustic Parameters	IEC	Units						
	$p_{ra,\alpha}$ at z_{MI}	(MPa)	3.39					
	P	(mW)		#		7.82		7.01
	$P_{1\times1}$	(mW)		7.67		7.67		
	z_s	(cm)			#			
	z_b	(cm)					#	
	z_{MI}	(cm)	0.60					
	$z_{pii,\alpha}$	(cm)	0.70					
	f_{awf}	(MHz)	5.19	5.09		5.09		5.19
Other Information	prr	(Hz)	#					
	srr	(Hz)	65.76					
	n_{pps}	-	11.00					
	$I_{pa,\alpha}$ at $z_{pii,\alpha}$	(W/cm ²)	493.59					
	$I_{spta,\alpha}$ at $z_{pii,\alpha}$ or $z_{sii,\alpha}$	(mW/cm ²)	76.75					
	I_{spta} at z_{pii} or z_{sii}	(mW/cm ²)	102.24					
	p_r at z_{pii}	(MPa)	3.27					
Operating Control Conditions	TX Level	(dB)	0.00	0.00		0.00		0.00
	Focus	(mm)	10.00	20.00		20.00		10.00
a	This Index is not relevant to this operating mode.							
b	This transducer is not intended for transcranial or neonatal cephalic uses.							
c	This formulation for TIS is less than that for an alternate formulation in this mode.							
#	No data is provided for this operation condition since the maximum index value is not reported for the reason listed.							

Acoustic Output Reporting Table — Track 3, FDA 510(k) and IEC 60601-2-37

(Per transducer/mode that exceeds MI or TI value of 1.0)

Transducer model: **9VE4**Operating mode: **2D-mode**

Index Label			MI	TIS		TIB		TIC
				At surface	Below surface	At surface	Below surface	
Maximum Index Value			1.50	0.45		0.45		0.69
Index Component Value				0.45	0.45	0.45	0.45	
Acoustic Parameters	IEC	Units						
	$p_{ra,\alpha}$ at z_{MI}	(MPa)	2.90					
	P	(mW)		#		18.00		38.25
	P_{1x1}	(mW)		17.73		17.73		
	z_s	(cm)			#			
	z_b	(cm)					#	
	z_{MI}	(cm)	1.10					
	$z_{pii,\alpha}$	(cm)	1.10					
	f_{awf}	(MHz)	3.82	5.34		5.34		4.89
Other Information	p_{rr}	(Hz)	#					
	s_{rr}	(Hz)	44.37					
	n_{pps}	-	2.00					
	$I_{pa,\alpha}$ at $z_{pii,\alpha}$	(W/cm ²)	362.93					
	$I_{spta,\alpha}$ at $z_{pii,\alpha}$ or $z_{sii,\alpha}$	(mW/cm ²)	13.87					
	I_{spta} at z_{pii} or z_{sii}	(mW/cm ²)	19.54					
	p_r at z_{pii}	(MPa)	3.23					
Operating Control Conditions	TX Level	(dB)	0.00	0.00		0.00		0.00
	Focus	(mm)	15.00	10.00		10.00		80.00
a This Index is not relevant to this operating mode.								
b This transducer is not intended for transcranial or neonatal cephalic uses.								
c This formulation for TIS is less than that for an alternate formulation in this mode.								
# No data is provided for this operation condition since the maximum index value is not reported for the reason listed.								

Acoustic Output Reporting Table — Track 3, FDA 510(k) and IEC 60601-2-37

(Per transducer/mode that exceeds MI or TI value of 1.0)

Transducer model: **9VE4**Operating mode: **M-mode**

Index Label			MI	TIS		TIB		TIC
				At surface	Below surface	At surface	Below surface	
Maximum Index Value			1.50	0.51		1.13		0.59
Index Component Value				0.51	0.32	0.27	1.13	
Acoustic Parameters	IEC	Units						
	$p_{ra,\alpha}$ at z_{MI}	(MPa)	3.00					
	P	(mW)		#		14.16		21.92
	$P_{1\times1}$	(mW)		21.92		14.16		
	z_s	(cm)			1.40			
	z_b	(cm)					1.30	
	z_{MI}	(cm)	1.50					
	$z_{pII,\alpha}$	(cm)	1.50					
	f_{awf}	(MHz)	4.06	4.89		3.94		4.89
Other Information	p_{rr}	(Hz)	720.05					
	s_{rr}	(Hz)	#					
	n_{pps}	-	#					
	$I_{pa,\alpha}$ at $z_{pII,\alpha}$	(W/cm ²)	571.76					
	$I_{spta,\alpha}$ at $z_{pII,\alpha}$ or $z_{sII,\alpha}$	(mW/cm ²)	323.00					
	I_{spta} at z_{pII} or z_{sII}	(mW/cm ²)	480.31					
	p_r at z_{pII}	(MPa)	3.58					
Operating Control Conditions	TX Level	(dB)	0.00	0.00		0.00		0.00
	Focus	(mm)	20.00	80.00		20.00		80.00
a	This Index is not relevant to this operating mode.							
b	This transducer is not intended for transcranial or neonatal cephalic uses.							
c	This formulation for TIS is less than that for an alternate formulation in this mode.							
#	No data is provided for this operation condition since the maximum index value is not reported for the reason listed.							

Acoustic Output Reporting Table — Track 3, FDA 510(k) and IEC 60601-2-37

(Per transducer/mode that exceeds MI or TI value of 1.0)

Transducer model: **9VE4**Operating mode: **Pulsed Wave Doppler**

Index Label			MI	TIS		TIB		TIC
				At surface	Below surface	At surface	Below surface	
Maximum Index Value			1.50	0.52		1.40		1.10
Index Component Value				0.52	0.33	0.43	1.40	
Acoustic Parameters	IEC	Units						
	$p_{ra,\alpha}$ at z_{MI}	(MPa)	3.22					
	P	(mW)		#		19.07		14.17
	P_{1x1}	(mW)		21.23		19.07		
	z_s	(cm)			1.30			
	z_b	(cm)					1.20	
	z_{MI}	(cm)	1.40					
	$z_{pii,\alpha}$	(cm)	1.40					
	f_{awf}	(MHz)	4.67	5.16		4.70		4.77
Other Information	p_{rr}	(Hz)	151.99					
	s_{rr}	(Hz)	#					
	n_{pps}	-	#					
	$I_{pa,\alpha}$ at $z_{pii,\alpha}$	(W/cm ²)	385.42					
	$I_{spta,\alpha}$ at $z_{pii,\alpha}$ or $z_{sii,\alpha}$	(mW/cm ²)	48.22					
	I_{spta} at z_{pii} or z_{sii}	(mW/cm ²)	75.51					
	p_r at z_{pii}	(MPa)	3.71					
Operating Control Conditions	TX Level	(dB)	0.00	0.00		0.00		0.00
	Focus	(mm)	80.00	80.00		20.00		5.00
a This Index is not relevant to this operating mode.								
b This transducer is not intended for transcranial or neonatal cephalic uses.								
c This formulation for TIS is less than that for an alternate formulation in this mode.								
# No data is provided for this operation condition since the maximum index value is not reported for the reason listed.								

Acoustic Output Reporting Table — Track 3, FDA 510(k) and IEC 60601-2-37

(Per transducer/mode that exceeds MI or TI value of 1.0)

Transducer model: **9VE4**Operating mode: **Color Flow Doppler**

Index Label			MI	TIS		TIB		TIC
				At surface	Below surface	At surface	Below surface	
Maximum Index Value			1.50	0.47		0.47		0.74
Index Component Value				0.45	0.45	0.45	0.45	
Acoustic Parameters	IEC	Units						
	$p_{ra,\alpha}$ at z_{MI}	(MPa)	3.35					
	P	(mW)		#		21.66		38.63
	P_{1x1}	(mW)		18.32		18.32		
	z_s	(cm)			#			
	z_b	(cm)					#	
	z_{MI}	(cm)	1.30					
	$z_{pii,\alpha}$	(cm)	1.30					
	f_{awf}	(MHz)	5.10	5.19		5.19		4.71
Other Information	prf	(Hz)	#					
	srr	(Hz)	2.33					
	n_{pps}	-	1.00					
	$I_{pa,\alpha}$ at $z_{pii,\alpha}$	(W/cm ²)	530.89					
	$I_{spta,\alpha}$ at $z_{pii,\alpha}$ or $z_{sii,\alpha}$	(mW/cm ²)	1.52					
	I_{spta} at z_{pii} or z_{sii}	(mW/cm ²)	2.32					
	p_r at z_{pii}	(MPa)	4.09					
	Operating Control Conditions	TX Level	(dB)	0.00	0.00		0.00	
Focus		(mm)	20.00	60.00		60.00		15.00
a This Index is not relevant to this operating mode.								
b This transducer is not intended for transcranial or neonatal cephalic uses.								
c This formulation for TIS is less than that for an alternate formulation in this mode.								
# No data is provided for this operation condition since the maximum index value is not reported for the reason listed.								

Acoustic Output Reporting Table — Track 3, FDA 510(k) and IEC 60601-2-37

(Per transducer/mode that exceeds MI or TI value of 1.0)

Transducer model: **CW2**Operating mode: **Continuous Wave Doppler**

Index Label			MI	TIS		TIB		TIC
				At surface	Below surface	At surface	Below surface	
Maximum Index Value			0.07	1.01		3.42		2.40
Index Component Value				1.01	0.80	1.01	3.42	
Acoustic Parameters	IEC	Units						
	$p_{ra,\alpha}$ at z_{MI}	(MPa)	0.10					
	P	(mW)		#		105.75		105.75
	P_{1x1}	(mW)		105.75		105.75		
	z_s	(cm)			1.70			
	z_b	(cm)					1.80	
	z_{MI}	(cm)	1.80					
	$z_{pii,\alpha}$	(cm)	1.80					
	f_{awf}	(MHz)	2.00	2.00		2.00		2.00
Other Information	prf	(Hz)	2000000.00					
	srr	(Hz)	#					
	n_{pps}	-	#					
	$I_{pa,\alpha}$ at $z_{pii,\alpha}$	(W/cm ²)	0.34					
	$I_{spta,\alpha}$ at $z_{pii,\alpha}$ or $z_{sii,\alpha}$	(mW/cm ²)	343.40					
	I_{spta} at z_{pii} or z_{sii}	(mW/cm ²)	541.47					
	p_r at z_{pii}	(MPa)	0.13					
Operating Control Conditions	TX Level	(dB)	0.00	0.00		0.00		0.00
	Focus	(mm)	#	#		#		#
a	This Index is not relevant to this operating mode.							
b	This transducer is not intended for transcranial or neonatal cephalic uses.							
c	This formulation for TIS is less than that for an alternate formulation in this mode.							
#	No data is provided for this operation condition since the maximum index value is not reported for the reason listed.							

Acoustic Output Reporting Table — Track 3, FDA 510(k) and IEC 60601-2-37

(Per transducer/mode that exceeds MI or TI value of 1.0)

Transducer model: **CW5**Operating mode: **Continuous Wave Doppler**

Index Label			MI	TIS		TIB		TIC
				At surface	Below surface	At surface	Below surface	
Maximum Index Value			0.05	1.07		2.38		1.24
Index Component Value				1.07	0.66	1.07	2.38	
Acoustic Parameters	IEC	Units						
	$p_{ra,\alpha}$ at z_{MI}	(MPa)	0.12					
	P	(mW)		#		44.86		44.86
	$P_{1\times1}$	(mW)		44.86		44.86		
	z_s	(cm)			1.40			
	z_b	(cm)					1.40	
	z_{MI}	(cm)	1.70					
	$z_{pii,\alpha}$	(cm)	1.70					
	f_{awf}	(MHz)	5.00	5.00		5.00		5.00
Other Information	p_{rr}	(Hz)	5000000.00					
	s_{rr}	(Hz)	#					
	n_{pps}	-	#					
	$I_{pa,\alpha}$ at $z_{pii,\alpha}$	(W/cm ²)	0.50					
	$I_{spta,\alpha}$ at $z_{pii,\alpha}$ or $z_{sii,\alpha}$	(mW/cm ²)	501.99					
	I_{spta} at z_{pii} or z_{sii}	(mW/cm ²)	934.11					
	p_r at z_{pii}	(MPa)	0.17					
Operating Control Conditions	TX Level	(dB)	0.00	0.00		0.00		0.00
	Focus	(mm)	#	#		#		#
a	This Index is not relevant to this operating mode.							
b	This transducer is not intended for transcranial or neonatal cephalic uses.							
c	This formulation for TIS is less than that for an alternate formulation in this mode.							
#	No data is provided for this operation condition since the maximum index value is not reported for the reason listed.							

Acoustic Output Reporting Table — Track 3, FDA 510(k) and IEC 60601-2-37

(Per transducer/mode that exceeds MI or TI value of 1.0)

Transducer model: **CW8**Operating mode: **Continuous Wave Doppler**

Index Label			MI	TIS		TIB		TIC
				At surface	Below surface	At surface	Below surface	
Maximum Index Value			0.04	0.62		1.21		1.02
Index Component Value				0.62	0.45	0.62	1.21	
Acoustic Parameters	IEC	Units						
	$p_{ra,\alpha}$ at z_{MI}	(MPa)	0.10					
	P	(mW)		#		16.29		16.29
	P_{1x1}	(mW)		16.29		16.29		
	z_s	(cm)			0.60			
	z_b	(cm)					0.60	
	z_{MI}	(cm)	0.60					
	$z_{pii,\alpha}$	(cm)	0.60					
	f_{awf}	(MHz)	8.00	8.00		8.00		8.00
Other Information	p_{rr}	(Hz)	8000000.00					
	s_{rr}	(Hz)	#					
	n_{pps}	-	#					
	$I_{pa,\alpha}$ at $z_{pii,\alpha}$	(W/cm ²)	0.34					
	$I_{spta,\alpha}$ at $z_{pii,\alpha}$ or $z_{sii,\alpha}$	(mW/cm ²)	341.40					
	I_{spta} at z_{pii} or z_{sii}	(mW/cm ²)	475.45					
	p_r at z_{pii}	(MPa)	0.12					
Operating Control Conditions	TX Level	(dB)	0.00	0.00		0.00		0.00
	Focus	(mm)	#	#		#		#
a This Index is not relevant to this operating mode.								
b This transducer is not intended for transcranial or neonatal cephalic uses.								
c This formulation for TIS is less than that for an alternate formulation in this mode.								
# No data is provided for this operation condition since the maximum index value is not reported for the reason listed.								

ACUSON Juniper
ACUSON Juniper Select
Diagnostic Ultrasound System
Clinical Applications Guide

US



An updated User and Reference Manual may be available for your ultrasound system. Access the most recent manual from the Siemens Healthineers Document Library at doclib.siemens-healthineers.com.

Product Version 3.0

Software Version VC10

©2018-2025 Siemens Medical Solutions USA, Inc.

All Rights Reserved.

Date of first issue: 2025-08

All trademarks are the property of their respective owners.

Siemens Medical Solutions USA, Inc. reserves the right to change its products and services at any time. In addition, this publication is subject to change without notice.

Manufacturer

Siemens Medical Solutions USA, Inc.
Ultrasound
22010 S.E. 51st Street
Issaquah, WA 98029
U.S.A.
Phone: +1-888-826-9702
siemens-healthineers.com

Siemens Healthineers Headquarters

Siemens Healthineers AG
Siemensstr. 3
91301 Forchheim
Germany
Phone: +49 9191 18-0
siemens-healthineers.com

Contents

Chapter A1	Annotations Procedures for annotating an image with text, body markers, and arrows.
Chapter A2	Biopsy and Puncture Procedures Procedures of the biopsy (puncture) function, including needle path verification and procedures for the needle visualization function.
Chapter A3	Panoramic Imaging Procedures for scanning an extended field of view with the panoramic imaging feature.
Chapter A4	<i>four</i>Sight 4D Imaging and 3-Scape Imaging Explanation of the Advanced <i>four</i> Sight 4D technology feature and the 3-Scape real-time 3D imaging feature.
Chapter A5	Strain Imaging Procedures for displaying the relative stiffness of tissue as an elastogram during real-time 2D-mode imaging. Includes procedures for performing a measurement with the shadow function and for determining a strain ratio.
Chapter A6	Virtual Touch Applications Procedures and descriptions for using acoustic radiation force impulse (ARFI) technology and quantitative ultrasound (QUS) technology to obtain qualitative visualizations or quantitative measurements of tissue stiffness. Includes the following applications: ▪ Point Shear Wave Elastography (pSWE) ▪ Virtual Touch Ultrasound-Derived Fat Fraction (UDFF)
Chapter A7	Contrast Agent Imaging Procedures for activating and using the contrast agent imaging feature.
Chapter B1	Physiologic Function Procedures of the physiologic function, including the ECG function and a troubleshooting guide.
Chapter B2	Stress Echo Procedures of the Stress Echo feature.

Chapter C1	Measurements and Calculations Descriptions of the measurement function with general procedures for measurements and calculations and the report feature.
Chapter C2	Obstetrical Measurements and Calculations Description of the measurement function with procedures for measurements and calculations and the report feature for Obstetrical exams.
Chapter C3	Cardiac Measurements and Calculations Description of the measurement function with procedures for measurements and calculations and the report feature for Cardiac exams.
Chapter D1	<i>syngo</i> Arterial Health Package Procedures of the <i>syngo</i> Arterial Health Package.
Chapter D2	eSie Left Heart Description of the software program for performing left heart volume measurements using 2D-mode data.
Chapter D3	Transesophageal Transducer Recommendations for leakage current testing, electrical safety, exam considerations, and cleaning and care information. Includes a description and illustrations of the articulation controls, on-screen indicators, messages for the transesophageal transducer, and temperature controls and safeguards.
Chapter D4	<i>syngo</i> Velocity Vector Imaging Procedures for activating and using the clinical software program to analyze cardiac clips. Includes an explanation of how the program tracks and estimates tissue velocity and other motion and deformation parameters at selected points on a user-defined outline (trace) of a structure.

Note: Not all features and options described in this publication are available to all users. Please check with your Siemens Healthineers representative to determine the current availability of features and options.

See also: For a list of available options and features, refer to Appendix A in the Instructions for Use.

About the User and Reference Manuals

The user and reference manuals contain descriptions for the following ultrasound systems:

- ACUSON Juniper diagnostic ultrasound system
- ACUSON Juniper Select diagnostic ultrasound system

Features and options unique to an ultrasound system are identified within the following locations in the user and reference manuals:

- Chapter 1, Instructions for Use
- Chapter 3, Instructions for Use
- Appendix A, Instructions for Use

The user and reference manuals consist of the following publications.

Publication	Includes
<i>Instructions for Use</i>	▪ Conventions and typographical conventions used in the manuals ▪ Intended audience ▪ Introduction to the ultrasound system ▪ Technical description of the ultrasound system ▪ Safety and care information for the system and compatible transducers ▪ Descriptions of system controls ▪ Procedures for system setup, examination fundamentals, image acquisition and optimization, and examination review ▪ Information about system security and protection of patient data against security threats ▪ Acoustic output data
<i>Clinical Applications Guide</i>	▪ Explanation of clinical applications and clinical software programs for use on the ultrasound system ▪ Procedures for the physiologic function and the biopsy function ▪ Description of the transesophageal transducer, including cleaning and care of the transducer ▪ General and exam-specific measurements and calculations ▪ Exam-specific annotations
<i>System Reference*</i>	▪ General and exam-specific measurement and calculation labels ▪ Description of customizable system settings ▪ Information about DICOM connectivity and network capabilities ▪ Clinical references ▪ Information regarding the electromagnetic compatibility (EMC) testing of this system

*Languages supported by the user interface include a translation of this publication.

Conventions

Take a moment to familiarize yourself with the following conventions in your user manual.

Convention	Description
Warnings and Cautions	 WARNING: Warnings are intended to alert you to the importance of following the correct operating procedures where risk of injury to the patient or system user exists.  Caution: Cautions are intended to alert you to the importance of following the correct operating procedures to prevent the risk of damage to the ultrasound system.
Notes and Prerequisites	Prerequisite: Prerequisites contain tasks the user must complete or information the user needs prior to performing a procedure. Note: Notes contain information concerning the proper use of the ultrasound system or correct execution of a procedure.
Examples of Cross-References	Examples: See also: For information about configuring a network, refer to Chapter 2 in the System Reference. See also: Touch screen controls, page 1-17 See also: For a list of transducers compatible with the needle guide accessories, refer to Chapter 6 in this manual.

Typographical Conventions

The following typographical conventions are in descriptions and procedures within your user and reference manuals. Use the conventions to identify the location of a control or selection on the ultrasound system.

Control or Selection	Description
Controls on the Control Panel and Keys on the Keyboard	Controls on the control panel are indicated by uppercase, boldface type. The function of the trackball and trackball keys is indicated by on-screen labels displayed on the image screen. Examples: ▪ Push DEPTH / FOCUS. ▪ Rotate 2D to adjust the 2D gain. ▪ Press CALIPER to activate the measurement function. ▪ Press and hold the Alt key on the keyboard to access the special characters. ▪ Press UPDATE. ▪ Position the marker and then press SET.
Controls on the Control Panel and Soft Key Selections on the Touch Screen	The touch screen has unlabeled controls located on the control panel. Tap, press, rotate or press and then rotate the control to make a selection or change a setting. The selection assigned to these controls is indicated by boldface type and enclosed in square brackets ([]). Examples: ▪ Press [Tx Power] downward to decrease transmit power. ▪ Rotate [Map].
Selections on the Image Screen	On-screen selections are indicated by boldface type. Roll the trackball to position the pointer on the on-screen control or object and then press a trackball key. Examples: ▪ Click Cancel. ▪ Double-click a thumbnail. ▪ Select the required check box to search for scheduled procedures.
Selections on the Touch Screen	On-screen selections are indicated by boldface type. To access the selection, tap the touch screen. To change the direction of an object displayed on the touch screen, drag the object on the touch screen. Examples: ▪ Tap Stress Echo. ▪ Drag the blue dot to indicate the direction of the arrow. ▪ Swipe the page to select the next page.

A1 Annotations

Annotating with Text	2
Annotating with Directional Arrows	4
Annotating with Body Markers	5

Annotating with Text

You can annotate an image by using the keyboard to enter text or by selecting system-defined or user-defined text. You can annotate an image with text during imaging or during review.

Use the configuration settings to customize text for the annotation function for each exam type:

- Create user-defined text
- Customize the text libraries
- Enable auto complete. With auto complete, the system displays predefined text matching the characters you enter on the image screen.
- Delete text from the image screen on unfreeze
- Configure user-defined controls to access the libraries of modifiers. The libraries include additional descriptive text to use in conjunction with the primary text describing the anatomy.

To annotate an image with text without activating the annotation function:

- Use the keyboard to enter the required text.

Note: To position the text at a location other than the default location, use the directional arrow keys to position the cursor, and then use the keyboard to enter the text.

To activate the text function:

- Press **TEXT BODY MARKER**.

To	Do This
Define the default position for text on the image screen	1. Use the trackball to position the cursor on the image screen. Or, use the directional arrows on the keyboard to position the cursor. 2. Press [Set Home].
Position the cursor at the default position	▪ Press [Home].
Insert annotation text on an image using auto complete	1. Position the cursor on the image screen and then enter the text. The system displays predefined text matching the entered characters, if specified in the configuration settings. When the entered text matches predefined text, the text is highlighted in green. 2. To accept the predefined text, press enter on the keyboard. Note: If you reposition the cursor before accepting the predefined text, the system displays only the characters you entered.

To	Do This
Insert predefined or user-defined annotation text on an image	Note: Each exam type has libraries of text for anatomical structures, right and left sides, and acquisition planes. 1. Position the cursor on the image screen. 2. Tap the selection for the required text.
Reposition text on the image screen	▪ Position the cursor on the text and then press SET. The system highlights the text in green.
Select text on the image screen	▪ Position the cursor on the text and then press SET.
Insert or remove characters in text on the image screen	1. Position the cursor between the characters. 2. Use the keyboard to insert or delete characters. The system displays editable text in white.
Delete one word of text on the image screen	1. Double-click the word. 2. Press [Delete Word].
Delete a line of text on the image screen	1. Click anywhere on the line of text. 2. Press [Delete Line].
Temporarily hide annotations	▪ To hide annotations, rotate [Hide Text].
Change the font size of the text	▪ Rotate [Font Size] and then select the required setting.
Remove all text, arrows, and body markers from the image screen	▪ Tap Clear Screen.
Exit the annotation function	▪ Press TEXT BODY MARKER.

Annotating with Directional Arrows

You can annotate an image with arrows during imaging or during review.

To activate the arrow function:

1. Press **TEXT BODY MARKER**.
2. Tap **Arrow**.

To	Do This
Select an arrow	▪ Select an arrow including type, color, and size in Arrow using the keys on the touch screen.
Change the angle of the arrow	▪ Drag the blue dot to indicate the direction of the arrow. ▪ Rotate ANGLE. ▪ Rotate [Rotate]. The system adjusts the angle of the arrow in 30° increments.
Insert an arrow	1. Position the cursor on the image screen. 2. Tap the Arrow tab. 3. Press SET.
Delete arrows	▪ Press [Delete].
Remove all text, arrows, and body markers from the image screen	▪ Tap Clear Screen.
Exit the arrow annotation function	▪ Press TEXT BODY MARKER.

Annotating with Body Markers

A body marker is an anatomical pictogram of anatomy in evaluation. You can add a transducer marker to indicate orientation of the transducer.

You can annotate an image with a body marker during imaging or review.

Note: The system retains the location of a body marker added during imaging when you change to a dual or quadrant format.

Use the configuration settings to customize body markers for each exam type.

To activate the body marker function:

1. Press **TEXT BODY MARKER**.
2. Tap **Body Marker**.

The system displays the first body marker on the annotation menu for the active exam type. The body marker includes a transducer orientation indicator.

To	Do This
Select a different body marker	▪ Tap the required body marker. Or, rotate SELECT. The body marker displays on the image screen and the touch screen.
Change the size of the body marker and transducer marker	(Available only for the active image) ▪ Rotate [Marker Size].
Reposition the body marker and transducer marker	1. Press [Trackball] to select Body Marker, if necessary. 2. Roll the trackball to reposition the markers.
Reposition only the transducer marker	▪ Drag the transducer marker on the touch screen to the required position. ▪ Tap the required position to locate the transducer marker. ▪ Use the trackball to position the transducer marker. a. Press [Trackball] to select Transducer Marker. b. Press SET to select the required option. ▪ Position indicates the trackball is assigned to repositioning the transducer marker. ▪ Rotate indicates the trackball is assigned to rotating the transducer marker. c. Roll the trackball to position the transducer marker.
Rotate the transducer marker	▪ To adjust the angle of the transducer marker in 15° increments, select a method: – Drag the blue dot on the touch screen to indicate the direction of the transducer marker. – Assign the trackball to rotate the transducer marker and then roll the trackball. – Rotate ANGLE. ▪ To adjust the angle of the transducer marker in 90° increments, press ANGLE.

To	Do This
Define the default position for the body marker on the image screen	1. Use the trackball to position the body marker on the image screen. 2. Press [Set Home].
Position the body marker at the default position	▪ Press [Home].
Delete a body marker from the image screen	▪ Press [Delete].
Remove all text, arrows, and body markers from the image screen	▪ Tap [Clear Screen].
Exit the body marker function	▪ Press TEXT BODY MARKER.

A2 Biopsy and Puncture Procedures

Biopsy (Puncture) Guideline Function	2
Activating the On-screen Guidelines	3
System Biopsy Safeguards.....	4
Needle Path Verification.....	5
Enhancing Visualization of a Needle	6

Biopsy (Puncture) Guideline Function

Refer to the in-box instructions for attachment and care procedures for all needle guide accessories, including information on reprocessing the reusable bracket.

⚠ WARNING: Percutaneous procedures always involve heightened risk to the patient and to the operator handling biopsy needle guides. Clinicians using Siemens Healthineers recommended biopsy devices under ultrasound guidance should be trained and must observe proper needle insertion sequencing with the needle guide in order to avoid undue discomfort and unnecessary risk and injury to the patient.

⚠ WARNING: The system displays biopsy guidelines to indicate an anticipated needle path. The guidelines do not account for a bent or deflected needle and are not intended as an absolute reference. It is the user's responsibility to verify correct positioning of the needle during a biopsy or puncture procedure.

⚠ Caution: Visually inspect the biopsy device before use. Do not attach a damaged biopsy device to a transducer.

Note: Linear and curved array transducers include reference lines as a supplementary visual guide during a freehand biopsy procedure. The reference lines are located at the center of the transducer surface and at points from the center.

The biopsy function is available for transducers compatible with needle guide attachments. The ultrasound system displays guidelines on the image to indicate an anticipated needle path. You must verify the angle of the needle guide on the biopsy device matches the angle selection on the image screen before performing any procedure.

You can use touch screen selections to display on-screen guidelines for use with transducers compatible with needle guide attachments. Use the configuration settings to assign activation of a biopsy or puncture procedure to a user-defined control on the control panel.

Activating the On-screen Guidelines

The biopsy function is activated only during real-time imaging in the following modes:

- 2D-mode full field of view (FOV)
- 2D-mode with Color flow Doppler
- 2D-mode with contrast

To activate the Guideline function:

Prerequisite: Refer to the needle guide manufacturer's instructions for attachment and care procedures for all needle guide accessories.

Note: Before performing any patient procedure using a needle guide, you must verify the needle path.

1. Tap **Biopsy** and then tap the required angle setting for the needle guide attached to the transducer.

The system displays an on-screen message as a precaution.

Message	Displays When
Verify that the physical needle guide matches your angle selection.	Full field of view (FOV) is active. Verify the angle selected on the needle guide attached to the transducer matches the angle selection on the touch screen.
Image has been changed to full FOV for Biopsy. Verify that the physical needle guide matches your angle selection.	Full field of view (FOV) was not active. The ultrasound system changed the image size to full field of view prior to activating the feature. Verify the angle selected on the needle guide attached to the transducer matches the angle selection on the touch screen.

2. Ensure the selected angle setting on the touch screen corresponds to the needle guide attached to the transducer and then click **OK**.
3. To change the angle setting, tap **Biopsy** and then tap a different setting.
Available angle settings are transducer dependent.
4. To add another guideline on the image screen, tap **Biopsy Line**.
5. Confirm the target of interest is indicated by the guidelines on the image screen.
6. To exit the biopsy function and remove the guidelines from the image, tap **Biopsy** and then tap **Off**.

System Biopsy Safeguards

 **WARNING:** To reduce the risk of injury to the patient, conduct a biopsy procedure during real-time imaging. When an image is frozen, the correct positioning of the needle cannot be assured.

While using the biopsy function, you can switch to a different acceptable mode or freeze the image. When you freeze the image, the color of the guidelines changes from yellow to white.

If you request an action that is not allowed during the biopsy function, the system displays a message indicating the action is not allowed.

Disconnecting the active transducer while in biopsy causes the system to exit the biopsy function.

Needle Path Verification

Before performing any patient procedure using a needle guide, you must verify that the path of the needle is accurately indicated by the on-screen guidelines.

The needle guide is ready for patient use **only after the path of the needle has been verified**.

Checklist of Items Required for Needle Path Verification:

- Transducer with attached needle guide
- Water-based coupling agent (gel)
- Sterile transducer cover
- New, straight, biopsy needle
- Sterilized container of sterilized and degassed water

To verify the path of the needle:

 **WARNING:** The system displays biopsy guidelines to indicate an anticipated needle path. The guidelines do not account for a bent or deflected needle and are not intended as an absolute reference. It is the user's responsibility to verify correct positioning of the needle during a biopsy or puncture procedure.

 **WARNING:** Do not use a needle guide if the path of the needle is not accurately indicated by the on-screen guidelines. The path of the needle must display within the guideline. Contact your Siemens Healthineers service representative if the needle path is not accurately indicated.

1. Attach the needle guide to the transducer.
2. Connect the transducer to the system and activate the transducer.
3. Set the system to the depth of the intended puncture procedure.
4. Tap **Biopsy**.

The system displays the following message as a precaution:

Verify that the physical needle guide matches your angle selection.

5. Ensure the selected angle corresponds to the needle guide attached to the transducer and then click **OK**.
6. To change the needle guide angle, tap **Biopsy**.
7. Immerse the head of the transducer into the degassed water and insert the needle into the needle guide.
8. Verify that the path of the needle displays according to the guidelines shown on the image screen.

After verification, the needle guide is ready for use.

Enhancing Visualization of a Needle


⚠ WARNING: Percutaneous procedures always involve heightened risk to the patient and to the clinician handling the needle. Clinicians should be trained in ultrasound needle guidance techniques.

⚠ WARNING: To avoid undue discomfort and unnecessary risk and injury to the patient, observe correct positioning of the needle throughout the procedure. Ensure the needle is in the enhanced image area. If the needle is in the unenhanced image area, do not continue the needle procedure until correct needle positioning is achieved.

⚠ WARNING: To improve visualization and avoid misinterpretation of the image, ensure the needle is shown perpendicular to the direction of the steered acoustic beam in the enhanced image area. If the needle is not inserted perpendicular, then the needle image may be less enhanced.

During a biopsy or interventional procedure using a needle with the linear transducer, you can select an angle setting for steering the acoustic beam of the transducer to enhance the visualization of the needle image. Visualization of the needle is optimal when the needle is perpendicular to the acoustic beam.

Other factors influence visualization of the needle, including the transducer and imaging parameter settings, needle gauge, insertion angle, tissue type, and user technique.

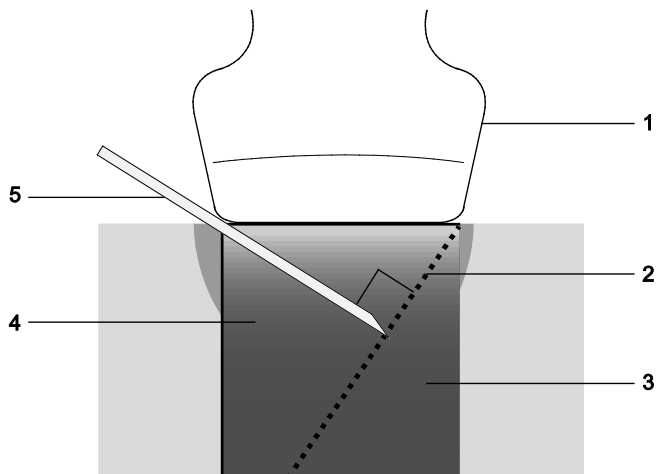

Example of needle visualization.

- 1 Linear transducer
- 2 Direction of the steered acoustic beam (depicted by a dotted line)
- 3 Unenhanced image area (no blue border)
- 4 Enhanced image area (depicted by a blue border)
- 5 Needle

To enhance visualization of a needle during real-time imaging:

(Not available for an image rotated in 90-degree increments)

Note: Use a recommended needle size of 16 to 20 gauge and thoroughly read and understand the warnings and descriptions beginning on page A2-6.

1. Tap **Needle Visual**. Or, press a user-defined control, if specified in the configuration settings.

Note: The system disables the following features and functions during needle visualization: compounding, changing to a trapezoid image format, and adjusting the angle of the field of view.

2. Select the intended trajectory of the needle.
 - **Angle Shallow** for a needle in a shallow depth.
 - **Angle Medium** for a needle in a mid-range depth.
 - **Angle Deep** for a needle in a deeper depth.

The system changes the angle of the dotted line on the image.

3. To select the enhancement level of the needle image, rotate **[Needle Enhance]**.
4. To reverse the horizontal (left/right) orientation of the direction of the steered acoustic beam, tap **L/R Flip**.

A3 Panoramic Imaging

Overview.....	2
Example Screen	2
Acquiring a Panoramic Image	4
Avoiding Artifacts During Acquisition	6
Measuring a Panoramic Image.....	7


Overview

⚠ WARNING: To avoid technique-related artifacts and measurement inaccuracies, read this entire chapter before using the Panoramic imaging feature.

Panoramic imaging acquires two-dimensional ultrasound images with a composite, extended field of view. Use CINE to review the image composed of hundreds of single frames.


Example Screen

Image creation begins with a standard field of view and extends the view during steady, forward motion scanning.

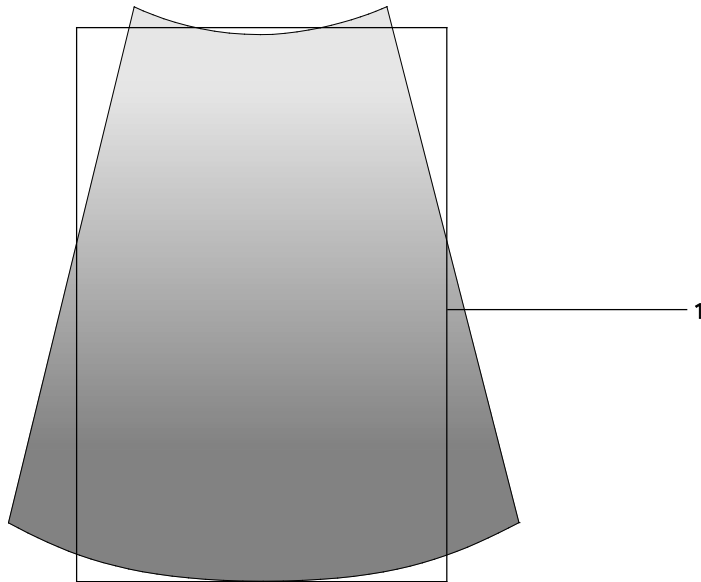

Example of an image screen during Panoramic setup.

- 1 Bounding box
Section of the scanned image used to create the composite image. If necessary, decrease the image depth to increase the image size so the image fills the bounding box.

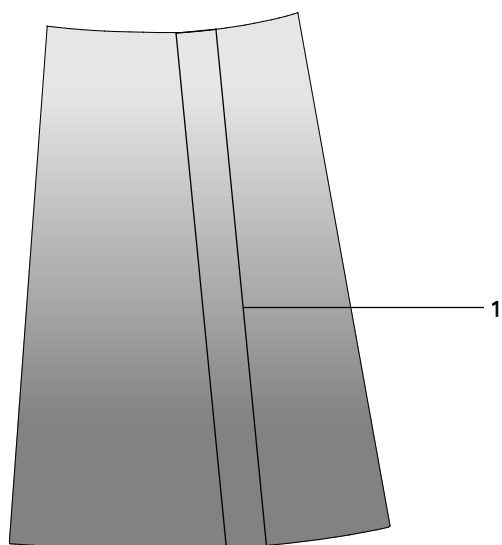

Example of an image screen during Panoramic image acquisition.

- 1 Region of Interest (ROI)
Boundary between the composite image and the real-time image.

Acquiring a Panoramic Image

Technique Tip: Perform a preview sweep of the scan plane in standard 2D-mode. The maximum length of a composite image depends on the selected depth. You can acquire approximately 8-10 times the selected depth for the length of the composite image for a linear scan.

To acquire a Panoramic image:

1. Optimize the 2D image.
Note: You cannot adjust the number of focal zones or persistence setting during image acquisition. The focal zone markers do not display on the composite image.
2. Position the single focal zone according to the transducer type:
 - Linear array transducers: lower half of image
 - Curved array transducers: upper half of image
3. Tap **Panoramic** and then press **SET** to begin acquisition.
 - **Acquire** indicates the trackball keys are assigned to beginning the acquisition.
 - **End** indicates the trackball keys are assigned to ending the acquisition.
4. Slowly and steadily scan in the required direction.
Note: As you scan, ensure the system displays the anatomy and structures as contiguous on the image screen.

 When the memory buffer for the composite image is full, the system stops acquisition and freezes the image.
5. To suspend or resume image acquisition, press **UPDATE**.
 - **Pause** indicates the control is assigned to suspending acquisition.
 - **Continue** indicates the control is assigned to resuming acquisition.

6. To stop image acquisition before the buffer is full, press **FREEZE**.

To	Do This
Review composite images in CINE review	▪ Tap Cine.
Display a flexible ruler with markers along the skin line during review or CINE review	▪ Tap Ruler. Or, use the configuration settings to enable the display of the ruler.
Resize a composite image after acquisition	▪ To scale the image to full acquisition size, tap Full Fit. ▪ To scale the image to fit the image area, tap Best Fit. ▪ To adjust the magnification of the image, rotate ZOOM. ▪ To scale the image in the size and rotation selected prior to entering CINE, tap Redisplay.
Pan or rotate a composite image	1. Press SET. – Pan indicates the trackball is assigned to panning the image. – Rotate indicates the trackball is assigned to rotating the image. 2. Roll the trackball. Note: You cannot pan an image when Best Fit is selected.
Rotate a frozen, composite image	▪ Rotate ANGLE clockwise or counterclockwise. Note: If Best Fit is selected, rotating the image causes the system to adjust the size of the image.
Activate standard 2D-mode	▪ Press 2D.
Begin acquisition of a new, composite image	▪ Tap Restart.
Set the composite image to the initial mode	▪ Tap Reset to set the frozen, rotated, full-fit, or best-fit composite images back to the default display mode.
Redisplay the previously acquired image	▪ Tap Redisplay.

Avoiding Artifacts During Acquisition

Artifact	Technique Tips
- Contour of the skin line has a "waving flag" appearance with folds - Jagged, broken skin line - Abrupt changes in boundaries or boundaries are not aligned - A structure visible in the real-time image abruptly disappears at the composite boundary	The transducer has moved off plane or off the required path. Stop your forward scanning motion and correct your orientation in the real-time 2D image before continuing acquisition of the composite image. To confine the image to a single plane or correctly following a vessel, observe the trailing edge of the real-time image area, located in the region of interest. Use the remainder of the real-time image as a guide to correct the rotation of the transducer while scanning in a forward motion.
Shadowing and bending	Balance the gain throughout the image to avoid a bending artifact. Low gain in the far field decreases the amount of data needed to ensure a geometrically correct composite image. Keep the transducer perpendicular and in full contact with the skin surface to avoid excessive bending and ensure a constant elevation position. A forward scanning motion avoids an image-compounding effect, which appears as a brighter, distorted section of the composite image. Scanning over a long length of superficial bone, occupying the full width of a single frame, creates shadowing. There is not enough data in the lower half of the image, causing excessive bending of the composite image. Scan a flat or gently curved surface for optimal results.
Swirl (image overlaps in depth)	Avoid performing a tightly curved scan (for example, transverse slice through the calf) may cause a swirl compounding artifact where the image overlaps in depth. Set the depth as shallow as possible. The depth should be less than the radius of the target being scanned. This leaves a blank "hole" in the image. If the depth is too deep, the bottom of the image appears swirled. Do not measure an image with a swirl at the bottom. Optimal Good Bad  Large Hole "Hole" Swirl Artifact <i>Example of curved Panoramic images.</i>

Measuring a Panoramic Image

 **WARNING:** To ensure accuracy, measurements should only be made on Panoramic images acquired in a single scan plane.

You can make general 2D-mode measurements on a frozen, full-sized, or best fit composite image. If the structure requiring measurement is contained within the boundary of a single, standard 2D-mode frame, make the measurement on the corresponding frame retrieved from CINE and not on the Panoramic image.

Technique Tips: Avoid measurements on the following composite images:

- Images not following a single plane, for example, when tracking a tortuous vessel. An out-of-plane image depicts the contour pattern of the skin line and the appearance of internal structures.
- Images with a swirl artifact at the bottom
- Across a large shadow in a Panoramic image
- Across gaps in an image, such as those encountered with a tightly curved scan
- Across a disrupted area of the image, indicated by a jagged flexible ruler along the skin line

Note: The system configuration option to activate the measurement function with the **FREEZE** key is not available during Panoramic imaging.

See also: For information on the measurement function, refer to Chapter C1 in this manual.

A4 *four*Sight 4D Imaging and 3-Scape Imaging

Overview.....	2
Activating 3D or 4D.....	2
Acquiring Volumes.....	3
Reviewing Volumes Stored in the CINE Buffer.....	4
Viewing Volumes and MultiPlanar Reformatting	5
Changing the Quadrant Display.....	5
Defining Volume of Interest and View Planes.....	7
Changing the Orientation of Volumes and MPRs	9
Visualizing the Anatomy of Interest in a Volume.....	10
Repositioning the Light around the Volume.....	11
Displaying Parallel Slices.....	12
Improving Contrast with Thick Slice Imaging.....	12
Straightening Curved Anatomy with Curved MPR.....	13
Removing Data from a Volume.....	14
Creating Presets for Volume Imaging.....	16
Measuring 3D and 4D Images	17
Storing and Printing Volumes and Images.....	17
Exporting Volumes or Clips	18
Imaging Techniques	18

Overview

 **WARNING:** To avoid technique-related artifacts, you must be thoroughly familiar with the techniques of 3D and 4D imaging. Read this entire chapter before using the 3-Scape real-time 3D imaging or *fourSight* 4D ultrasound imaging technology.

Ultrasound data based on three-dimensional imaging methods can facilitate the diagnostic process. You can view and manipulate the volumes to differentiate anatomic structures (for example, by selecting an arbitrary slice, modifying volume rendering, or rotating and magnifying localized areas of interest).

The 3-Scape imaging feature allows the acquisition of three-dimensional ultrasound images.

The *fourSight* 4D imaging feature allows the acquisition and simultaneous viewing of three-dimensional ultrasound images in real-time. You can use the feature to access motion, such as fetal movement.

MultiPlanar Rendering (MPR) provides a view of each plane of the volume (acquisition, elevation, and coronal) as an arbitrary slice.

See also: For information on transducers compatible with 3D or 4D imaging, refer to Appendix A in the Instructions for Use manual.

Activating 3D or 4D

The ultrasound system activates 3D or 4D imaging depending on the selected transducer and exam type, or as specified in the configuration settings.

To activate 3D or 4D:

Prerequisite: Optimize the 2D-mode image before beginning 3D or 4D imaging acquisition.

1. Press **3D / 4D**.
2. Tap **3D** or **4D**.
3. To exit, press **3D / 4D** again. Or, press **2D**.

Acquiring Volumes

During 3D volume acquisition, the system updates the incremental data set or slice (acquired image). The incremental slice is the mid-plane perpendicular to the 2D scan plane. Use the incremental slice to assess imaging artifacts and the position within the selected scan angle. Use the time gauge to estimate an optimal scanning speed. The time gauge depicts the elapsed time, scan length or angle, and scan duration.

During 4D imaging, the system displays one volume at a time and displays the acquisition rate in volumes per second (vps).

Note: At any time during acquisition and viewing, you can pause volume acquisition to review the data sets stored in the CINE buffer.

To acquire a volume:

Prerequisite: Optimize the 2D-mode image before beginning 3D or 4D imaging acquisition.


1. Press **3D / 4D**.
2. Tap **3D** or **4D**.
3. Size and position the region of interest.
4. Select the acquisition angle and quality.
 - Rotate **[Angle]** to adjust the scan angle.
 - (Available only for 4D imaging)
 - Rotate **[Steering Angle]** to adjust the transducer scan angle.
 - Rotate **[Quality]** to adjust the volume acquisition quality.

To	Do This
Select a scan method	(Available only for 3D imaging) 1. Select the scan method: – Rocked requires an angular scanning motion. Use the transducer face as a pivot. – Linear requires a horizontal or vertical sweep of the transducer. 2. To adjust the scan speed, rotate [Speed] . The system determines the scan duration based on your selections and the current frame rate.
Begin acquisition	Note: Hold the transducer steady during acquisition. ▪ Press UPDATE. An audible tone indicates volume acquisition has begun. A second audible tone indicates volume acquisition is complete. The system will display the volume of interest (VOI) within the region of interest.

To	Do This
Temporarily pause and resume volume acquisition	▪ Press FREEZE .
Begin another volume acquisition without saving the displayed volume	1. Press FREEZE . 2. Press UPDATE .
Cancel the acquisition (for example, for an image artifact)	▪ Press UPDATE .
End 3D acquisition before the scan duration on the time gauge elapses	▪ Press FREEZE . For 3D imaging, the system processes the acquired information into one or more data sets and then displays the volume.
End 4D acquisition	▪ Press 3D / 4D .

Reviewing Volumes Stored in the CINE Buffer

When 4D volume acquisition is paused, you can review the volumes stored in the CINE buffer. The system displays a 4D CINE bar at the bottom of the screen when acquisition is paused.

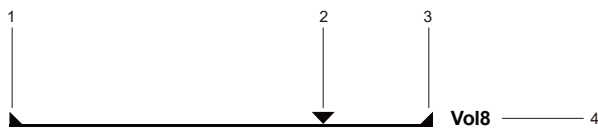

Example of 4D CINE bar.

- 1 Left CINE marker indicates the beginning of the CINE data.
- 2 Volume Indicator indicates the location of the volume within the CINE data and the direction of review (backward or forward).
- 3 Right CINE marker indicates the end of the CINE data.
- 4 Volume counter indicates the number of the displayed volume.

Viewing Volumes and MultiPlanar Reformatting

The system displays the volume and multiplanar reformatting during 4D volume acquisition, after 3D volume acquisition, or during review of retrieved volumes.

The MPRs are oriented to the center of the volume; each slice is orthogonal to the other two slices. The reference point, represented by a dot within each quadrant, indicates the point of intersection of the MPRs.

Changing the Quadrant Display

Quadrant D contains the volume. Quadrants A, B, and C contain MPRs of the volume; A is the acquisition plane, B is the elevation plane, and C is the coronal plane.


To	Do This
Select a quadrant	Position the pointer on the required quadrant and then press SET.
Cycle through quadrants	Tap the required quadrant.
Display or hide the volume quadrant	- To hide the volume, tap Render and then tap MPR Only. - To display the volume, tap Render and select a setting other than MPR Only.
Select a display method for the volume quadrant	- Tap the required method: Slice displays a one-voxel-thick slice in the three-dimensional context. The displayed slice (within the volume) corresponds to the most recently selected slice (MultiPlanar Reformatting or MPR) quadrant. Surface smooths image contours, creating a soft, sculptured appearance for highlighting surface features. Uses the mode-specific selections for opacity percentage and thresholds. MPR Only increases the rendering speed by displaying the MPRs (arbitrary slices) only. MinIP (Minimum Intensity Projection) limits the display to minimum-value pixels by selecting pixels with the lowest value along each projection of the volume. May be helpful in visualizing hypo-echoic lesions and large vascular structures. MaxIP (Maximum Intensity Projection) limits the display to maximum-value pixels by selecting pixels with the highest value along each projection of the volume. May be helpful in visualizing hyper-echoic structures and fetal skeletons. MeanIP (Mean Intensity Projection) limits the display to average-value pixels by selecting pixels with the arithmetic mean (average) value along each projection of the volume.

To	Do This
Change the display format for the volume and MPRs (for example, 1:1 or 4:1)	▪ Tap Render and then tap the required format. – 1:1 displays the selected quadrant only  – 2:1 displays two quadrants on the screen  – 1:3 display the volume on the left of the screen and the other three quadrants on the right  – 4:1 displays all four quadrants on the screen 
Display or hide the boundaries of the volume	(Available only when the volume of interest is disabled) ▪ Tap Wireframe.

Defining Volume of Interest and View Planes

You can define the data to display in the volume quadrant by enabling the VOI (volume of interest) and then resizing and repositioning the VOI in relation to the acquired volume data.

The system displays the VOI in the volume and the related view plane in each orthogonal MPR. The view plane is the side of the VOI represented by a solid or dotted line.

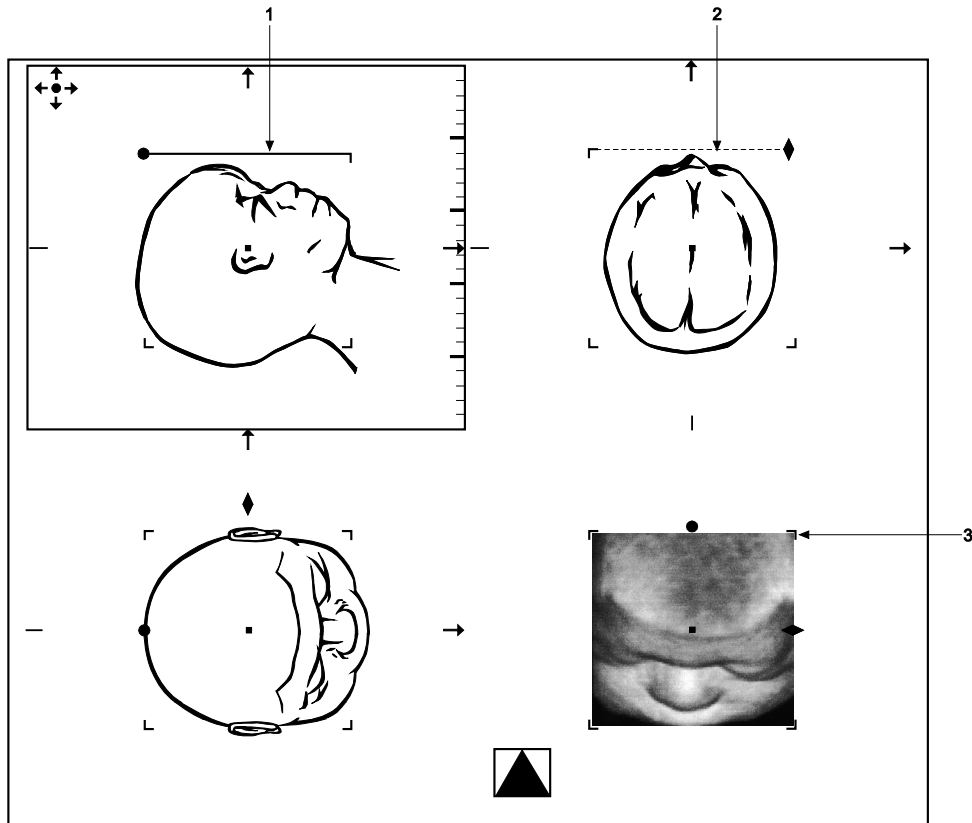

Example of the VOI and view planes.

- 1 View plane for the selected quadrant
- 2 View plane for an unselected quadrant
- 3 View of interest

To	Do This
Display or hide the VOI	Tap VOI and then select the required option: On, Off.
Change the size of the VOI	- Select a quadrant. - Press SET, if necessary. The system displays the VOI size icon.  - Roll the trackball to adjust the VOI size.
Reposition the volume	- Select the volume quadrant. - Tap VOI to activate the volume of interest, if necessary. - To activate the pan control, press SET. - Roll the trackball to pan the image. The system displays the pan volume icon. 
Reposition an MPR image	- Select a quadrant. - Press SET, if necessary. The system displays the pan icon.  - Roll the trackball to pan the image.
Curve the rendering view line of the VOI to exclude extraneous information	(Available only when an MPR is selected) - Tap Curved. - Roll the trackball to change the contours of the curve.
Change the render direction of the VOI	- Tap RenderDir. - Select the required direction to start the rendering: Up/Down begins rendering from the top Down/Up begins rendering from the bottom Front/Back begins rendering from the front Back/Front begins rendering from the back Left/Right begins rendering from the left Right/Left begins rendering from the right

Changing the Orientation of Volumes and MPRs

Rotating a volume or MPR along any of the three axes can help you visualize an anatomic structure. You can rotate a volume manually or automatically.

The reference orientation indicator in the volume quadrant uses color to indicate the relative location of the data acquired first (blue) and the data acquired last (red). The top (apex) of the pyramid represents the transducer. The base of the pyramid represents the data at the greatest depth.



Axis of rotation	Volume rotation icon	MPR rotation icon
X axis		
Y axis		
Z axis		

Note: The system displays a reverse symbol ("R") on the MPR quadrants with a reversed orientation in relation to the original acquisition.

To	Do This
Manually rotate the volume or MPR	- Press UPDATE to select the axis of rotation and then roll the trackball to rotate the volume or MPR. - Rotate the required control on the control panel: X rotates along the x axis Y rotates along the y axis Z rotates along the z axis The system displays a rotation icon on the selected MPR or volume. When an MPR is selected, the system also displays a line indicating the axis of rotation on the two corresponding MPRs.
Rotate the volume in 90° increments	Tap the required rotation angle, for example, 0, 90, 180, 270.
Enable automatic rotation of the volume	Note: You can change the automatic rotation settings while the volume is rotating. - Tap Animation. - Adjust the settings for speed, angle, axis, and clip time. - Tap the play icon. - To pause animation, tap the pause icon.

To	Do This
Restore the orientation from the original acquisition for all quadrants	▪ Tap Reset and then tap Reset Ori.
Restore the orientation and all imaging parameters from the original acquisition	▪ Tap Reset and then tap Reset All.
Reverse the render direction	▪ Tap Reset and then tap Flip.
Align the volume with the selected MPR	1. Select an MPR. 2. Tap Reset and then tap Sync.

Visualizing the Anatomy of Interest in a Volume

The cut plane function helps you visualize the anatomy of interest, such as a lesion, that is not visible on the surface of a volume by controlling the amount of data displayed in the volume.

For example, if anatomy is visible in an MPR but is not visible in the volume, you can use the cut plane function to reveal the anatomy of interest in the volume.

Note: For an optimal view of the cut plane, you can rotate the volume or align the orientation of the volume with that of the selected MPR.

To visualize the anatomy of interest in the volume:

Prerequisite: Disable the volume of interest to control the plane display in the volume.

1. Tap **CutPlane** and adjust the setting:
 - Select **A** to display the data in front of the cut plane.
 - Select **B** to display the data behind the cut plane.

The system removes all mode-specific volume data beyond the cut plane. A line on the reference orientation indicator indicates the location of the cut plane.

2. To activate the slicing function, roll the trackball.

Repositioning the Light around the Volume

Use the trackball or touch screen selections to rotate the direction of light 360° around the volume. The light illuminates surfaces and enhances shadows on the volume. You can adjust the light intensity to control the light reflection on the surface of the volume and to optimize the lighting effect.

To activate and position the direction of the light around the volume:

Prerequisite: Activate 3D or 4D imaging.

1. Tap **LightSource**.

The volume displays in a 1:1 format with an icon to indicate the feature is active.



2. To reposition the direction of light using the trackball:

- a. Press **UPDATE**.

A marker for the light source displays on the image screen.



- b. Roll the trackball to reposition the marker.

The system updates the direction of the light source on the volume.

3. To reposition the direction of light using the touch screen:

- a. Tap **Lighting>**.

- b. Tap the required touch screen selection.

Selection	Description
	Repositions the light source directly in front of the volume.
	Repositions the light source in the direction indicated by the arrow. Selections are available for multiple directions of light, including other angles around the volume.

4. To adjust the brightness of the light, rotate **[Light]**.

Displaying Parallel Slices

You can view multiple parallel slices of a plane at the same time. MultiSlice enables you to view anatomy of interest in sequential order, similar to a CT or MRI scan. You can also use MultiSlice to view a 4D CINE sequence.

To	Do This
Activate MultiSlice	▪ Tap MultiSlice .
Slice the volume in the selected MPR	▪ Tap Vertical or Horizontal . The volume quadrant displays a page of the slices. A blue dot indicates the center slice of the acquired volume.
View a slice in full screen format	▪ Press SELECT .
Display the next slice of the image	▪ Tap Next Image .
Display the previous slice of the image	▪ Tap Prev Image .
Display the next page of slices	▪ Tap Next Page .
Display the previous page of slices	▪ Tap Prev Page .
Change the division between pages	1. Select the quadrant. 2. Roll the trackball to scroll the images.
Adjust the distance between slices	▪ Rotate [Slice Spacing] .
Change the number of slices displayed on each page	▪ Tap the required layout.
Change the reference MPR (slice quadrant)	▪ Tap the required plane.

Improving Contrast with Thick Slice Imaging

Thick Slice Imaging (TSI) allows you to adjust the thickness of an individual slice to improve contrast resolution.

To	Do This
Activate TSI	▪ Tap TSI .
Display a thick slice image of the coronal plane (with the acquisition plane as reference)	▪ Tap A/C .
Display a thick slice image of the acquisition plane (with the acquisition plane as reference)	▪ Tap A/A .
Change the thickness of the TSI image	▪ Tap Slice Thickness and then select a setting.
Rotate the reference line on the reference plane around the z axis	(Available only for the A/C option) 1. Tap TSI . 2. Press SET to activate the rotation function. 3. Roll the trackball to the required position. The TSI image in the volume quadrant corresponds to the location of the reference line.
Exit TSI	▪ Tap Exit .

Straightening Curved Anatomy with Curved MPR

You can use curved MPR to straighten curved anatomy (such as a fetal spine) to facilitate evaluation.

Curved MPR generates cross-sections along a centerline to show details of a curved structure and displays the result with the orthogonal reference image of the displayed volume.

Note: Take care in drawing and interpreting the curved MPR result. Structures outside the path may be distorted because the final image is straightened. To avoid misrepresenting and distorting anatomical structures and their relationships, ensure that the path you draw does not cross the same object more than once.

To activate Curved MPR:

1. Optimize the orientation and display the curved anatomy in quadrant A.
2. Tap **Curved MPR**.

To	Do This
Define a line curve	1. Tap Line to draw a straight line.  2. Position and then anchor the starting point. 3. Position and then anchor the ending point.
Define a spline curve	1. Tap Spline to draw a line with curved segments.  2. Position and then anchor the starting point. 3. To insert individual points along the curve, roll the trackball and then press UPDATE at each interval. 4. Position and then anchor the ending point. The volume quadrant displays the curved anatomy as straightened. The arrows indicate the direction that the curve was drawn. A red dot indicates the center of the drawn line.
Trace a curve	1. Tap Trace.  2. Position and then anchor the starting point. 3. Roll the trackball to create a trace and then press SET. 4. Position and then anchor the ending point. The volume quadrant displays the curved anatomy as straightened. The arrows indicate the direction that the curve was drawn. A red dot indicates the center of the drawn line.
Reposition the drawn MPR curve to another location on the image	1. To display the repositioning pointer, position the cursor over the center segment of the line. The system displays a hand cursor. 2. Drag the line to the required location and then press SET. Or, press the arrow keys to reposition the line to the required location.
Rotate the MPR curve around the center point of the image	1. To display the rotation pointer, position the cursor over either endpoint of the line. 2. Drag the endpoint in the direction of the required rotation and then press SET.
Exit Curved MPR	▪ Tap Exit.

Removing Data from a Volume

You can remove data from a volume to clarify anatomical structures.

Note: To determine the type of edit required, rotate the volume.

To	Do This
Display the editing tools	▪ Tap Editing.
Remove data from inside or outside a polygon	Note: You can perform multiple polygon edits. 1. Tap Polygon.  The system constrains the cursor to the selected quadrant. 2. Position and then anchor the starting point. 3. For each segment of the polygon, roll the trackball to position the ending point and then press UPDATE. The system uses the ending point of the segment as the starting point for the next segment. 4. Press SET to complete the polygon. 5. Position the cursor where data will be removed (inside or outside the polygon) and then press SET.
Remove data between parallel planes	(Available only when the volume of interest is disabled) Note: You can simultaneously perform multiple parallel cut edits. 1. Tap Parallel Cut.  The system displays the wireframe and highlights the active plane. 2. For each plane to be used in the edit, roll the trackball to position the plane and then press UPDATE. 3. Press SET to complete the edit.
Remove data inside a niche (cube shape)	(Available only when the volume of interest is disabled) Note: You can perform one niche edit at a time. 1. Rotate the volume until the area that requires editing is in front. 2. Tap Niche.  The system outlines a cube in green on the volume. The cube represents the niche edit. 3. To simultaneously adjust the size of all planes of the cube, roll the trackball. 4. To resize one plane of the cube, press UPDATE. 5. Press SET to complete the edit.

To	Do This
Remove data inside or outside a trace	Note: You can perform multiple trace edits. 1. Tap Trace.  The system constrains the cursor to the selected quadrant. 2. Position and then anchor the starting point. 3. Roll the trackball to trace the area of interest. 4. Position and then anchor the ending point. 5. Position the cursor where data will be removed (inside or outside the traced shape) and then press SET.
Erase data in a volume or slice	Note: When you complete the edit, the system updates the other quadrants. 1. Select a tool for removing image data: – Tap Small Eraser for fine editing.  – Tap Large Eraser to remove large areas of background data.  The system displays a circle-shaped cursor on the most recently selected quadrant. 2. Position and then anchor the starting point. 3. Roll the trackball to erase data. 4. Press SET to complete the edit.
Complete the current edit	▪ Choose a method: – Select a different quadrant. – Magnify the selected quadrant. – Change the orientation of the selected quadrant.
Cancel the most recent segment of the current edit (such as an eraser edit)	▪ Rotate [Undo Last Edit].
Cancel all completed edits	Note: VOI edits cannot be canceled. ▪ Tap Undo All Edit.

Creating Presets for Volume Imaging

User-defined presets are an optimized configuration of imaging parameters for use during 3D or 4D imaging.

To create a preset for 3D or 4D:

1. Adjust the image optimization settings.
2. Tap **SubPreset** and then tap **Define**.
3. Click **New** and enter a name for the preset.
4. To enable the new preset as the default exam type for subsequent patient studies, select the **Set as default** check box.
The system displays an asterisk next to the name.
5. Click **OK**.
6. Click **Close**.
7. To begin using the new preset, click **Apply**.

To edit a preset for 3D or 4D:

1. Adjust the image optimization settings.
2. Tap **SubPreset** and then tap **Define**.

To	Do This
Enable an existing preset as a default exam type	1. Select a preset. 2. Click Default. The system displays an asterisk next to the name.
Rename an existing preset	1. Select a preset. 2. Click Rename and enter a new name. 3. Click OK.
Overwrite an existing preset	1. Select a preset. 2. Click Overwrite. 3. Click OK.
Delete a preset	Note: You cannot delete system-defined presets. 1. Select a preset. 2. Click Delete.

Measuring 3D and 4D Images

The basic measurement tools are valid for measurements on a quadrant.

Do not measure across the MPR quadrants or volumes. Perform each measurement in an individual MPR quadrant.

See also: For information about the measurement function, refer to Chapter C1 in this manual.

Storing and Printing Volumes and Images

Images are stored as single-frame images in screen capture format on the system. The images are stored within the series object that contains the volume. Series objects are accessible from the patient browser.

Note: An image of the volume can be printed. The volume itself (all of its multiple frames) cannot be printed.

Use the patient browser to transfer volumes (and their representative images) to connected devices that support multi-frame images.

To	Do This
Store a 3D or 4D volume	1. Press UPDATE. 2. Press the key assigned to the clip store function. 3. (Available only for 4D volumes) a. Press FREEZE. b. Select 4D Single volume and then click OK.
Store a 3D clip of the volume during automatic rotation	1. Tap Animation. 2. Tap the play icon. 3. Press the key assigned to the clip store function.
Temporarily store a 4D clip of the volumes in the CINE buffer	1. Press UPDATE to acquire a volume. 2. Press the key assigned to the clip store function. 3. Press FREEZE. 4. Select 4D Cine. 5. Click OK.

Exporting Volumes or Clips

You can export the volumes or clips of images to a storage location.

See also: For information about exporting a patient study, refer to Chapter 5 in the System Reference manual.

Imaging Techniques

Adequate Gel

Apply a generous amount of coupling agent (gel) along the entire area to be scanned to avoid disrupting the scan sweep. Insufficient gel causes the transducer to drag on the skin during a linear scan.

Preview Sweep

Before volume acquisition, perform a preview scan of the anatomy in standard 2D-mode. To be sure you acquire a complete structure of interest, include some of the surroundings in the 2D-mode image.

Scan Speed (for 3D imaging methods)

Scan at a constant rate. An irregular movement pattern may result in motion artifacts or distortions of the 3D image.

Transducer (Linear and Rocked Scan Methods)

To ensure a constant elevation position during a linear scan, hold the transducer so that your little finger is in contact with the skin alongside the transducer face. This acts as a guide to ensure that the transducer is parallel to and in full contact (perpendicular) with the skin surface.

For a rocked scan, ensure a steady, even scan motion.

A5 Strain Imaging

Overview.....	2
Activating Strain Imaging	2
Basic Measurements	3
Shadow Function.....	3
Strain Ratio.....	4
Calculated Results of the Strain Ratio Tool.....	5
Patents.....	5

Overview

Strain imaging is a qualitative visualization method used to display the relative stiffness of tissue as an elastogram.

The ultrasound system displays two images on the image screen in live dual format.

Activating Strain Imaging

To activate Strain imaging:

1. Tap **Virtual Touch** during 2D imaging.
The ultrasound system displays two images on the screen in a live, dual format. The left image is a 2D-mode image with a region of interest (ROI). The right image is a 2D image with a corresponding ROI showing the elastogram.

The system also displays a real-time quality factor (**QF**) score. A value greater than 60 indicates minimal global motion artifact. A value less than 50 indicates a less than optimal result; perform another acquisition.
2. To regenerate the elastogram, use cycles of gentle up and down compression while avoiding lateral and vertical movement of the transducer.
Note: The natural compressions generated by a patient's respirations may be sufficient to generate the elastogram.
3. To increase the size of the live dual format on the screen, rotate **ZOOM**.
As you enlarge the image, the system crops the edges of the image while maintaining the aspect ratio.
4. To view a full screen Strain image, tap **Live Dual**.
5. To apply a gray scale map, tap **EI Map Index** and then select an option.
 - **1** applies the standard map. This is the natural result of the strain estimation from the echo signal without compensation for pre-compression (stress) or non-uniform stress applied by the user.
 - **2** applies a normalized map. This compensates for the potential that non-uniform stress is applied with the transducer during scanning.
 - **3** applies an inverted map.

The system displays an elasticity imaging bar above the gray bar for the active map. The elasticity imaging bar represents the range of gray shades available for depicting the stiffness of the tissue from **SF** (soft) to **HD** (hard).
6. To apply a color map, tap **EI Color** and then select an option.
7. To exit, tap **Exit**.

Basic Measurements

You can perform all basic measurements including the area ratio and distance ratio during Strain imaging on either image in a dual screen layout.

See also: For information about the measurement function, refer to Chapter C1 in the Clinical Applications Guide.

Shadow Function

When the shadow function is active, the system duplicates and displays the measurement on the other image.

To perform a measurement using the shadow function:

1. Activate the measurement function.
2. Tap **Live Dual**, if necessary.
3. Press **CALIPER** to activate the measurement function.
4. Tap **Shadow**, if necessary.
5. Measure a distance, trace, or ellipse on either image.
6. To exit the shadow function, tap **Shadow**.

Strain Ratio

Strain imaging provides a compound measurement tool to compare the relative stiffness of tissue within two regions of interest, for example, compare the average strain of a lesion and the average strain of the surrounding tissue.

- A strain ratio of 1 indicates the stiffness of the tissue in both regions of interest is equal.
- A strain ratio greater than 1 indicates the tissue in the first region of interest (for example, the surrounding tissue) is softer than the tissue in the second region of interest (for example, a lesion).
- A strain ratio less than 1 indicates the tissue in the first region of interest is harder than the tissue in the second region of interest.

When interpreting the results of strain ratio calculations, consider the type of tissue within each region of interest. For example, results will be different depending on whether you selected the lesion or the surrounding tissue as the first region of interest.

You can label strain ratio measurements with any labels available with exams that support Strain imaging.

To determine the strain ratio:

Prerequisite: Select a linear transducer, activate Strain imaging, and select a map index of 1 or 2.

Note: The system calculates the strain ratio using only the region of interest displayed within the elastogram. If the strain ratio ROI extends beyond the elastogram, expand the elastogram. To expand the elastogram, press a trackball key to select **Size**, and then roll the trackball to adjust the size of the elastogram, and then press a trackball key again.

1. Activate the measurement function.
2. Tap **Strain Ratio**.

The system activates the strain ratio tool.

Note: If the strain ratio ROI extends beyond the elastogram, the system calculates the strain ratio using only the portion of the ROI within the elastogram. To prevent the ROI from extending beyond the elastogram, expand the elastogram to include the entire image.

3. Determine the first region of interest for reference.
 - a. Roll the trackball to position the first marker.
 - b. To resize the diameter of a ROI, rotate **[ROI Diameter]**.
 - c. Press **SET**.
4. Determine the second region of interest for a mass.
 - a. Roll the trackball to position the second marker.
 - b. To resize the diameter of a ROI, rotate **[ROI Diameter]**.
 - c. Press **SET**.

The system displays the strain ratio and the average strain (%) for both regions of interest in the measured results.

Calculated Results of the Strain Ratio Tool

Results	Abbreviation
Region of interest for reference (surrounding tissue)	ROI1 (Ref)
Region of interest for a mass (possible lesion)	ROI2 (Mass)
Strain ratio by comparing the relative stiffness of tissue within two user-selectable regions of interest	Strain Ratio

Patents

Patent numbers in the United States of America for the Strain imaging feature include the following numbers:

- 6,270,459

A6 Virtual Touch Applications

Acoustic Radiation Force Impulse Technology	2
Quantitative Ultrasound Technology	3
Virtual Touch Point Shear Wave Elastography (pSWE)	3
Virtual Touch Ultrasound-Derived Fat Fraction (UDFF)	4
Measurement Labels	5
Virtual Touch Patient Report and Left Panel	6
Scanning Techniques to Mitigate Technical Limitations	7
Shear Wave Velocity and Elasticity Measurements in Liver Assessment.....	9

Acoustic Radiation Force Impulse Technology

⚠ WARNING: To prevent the possibility of inadvertently heating tissue, avoid repeated image acquisitions near the interface of bone and tissue within a short period of time.

⚠ WARNING: Measurement values for shear wave velocity, elasticity, and other image display settings are relative indexes intended only for the purpose of comparison with other measurements performed using the ACUSON Juniper ultrasound system. Absolute values for shear wave velocity and elasticity measurements may vary among different manufacturers due to multiple system-dependent factors, including shear wave frequency, excitation beam (push beam) frequency, shear wave detection techniques, and shear wave velocity and elasticity estimation methods.

The system displays the following message on the screen during Point Shear Wave Elastography (pSWE) and in the report if shear wave speed measurements are saved to the report:

Note: Shear wave speed and elasticity values may vary among manufacturers!

The following features use acoustic radiation force impulse (ARFI) technology to obtain qualitative visualizations or quantitative measurements of tissue stiffness:

- Point Shear Wave Elastography (pSWE)

Note: The Ultrasound-Derived Fat Fraction feature does not use ARFI technology.

ARFI technology uses a time-controlled sequence of "push pulses" from the transducer, rather than manual compression or palpation cycles, to induce a small displacement of tissue. The image qualitatively depicts the relative stiffness of tissue within a selected region of interest (ROI) compared to surrounding tissue. The system can also provide a quantitative measurement of shear velocity (Vs) or elasticity (E) at a specific anatomical location.

During acquisition, ARFI technology combines a period of pulse sequences with a temporary freeze period of no transmit power. A timer indicates the remaining time in seconds until the transducer is available for imaging. This duration is typically a few seconds depending on the transducer and imaging conditions. The cooling timer prevents possible heating of the transducer surface from the momentary shear wave generation pulses.

Cooling 5

Example of Virtual Touch timer.

The system uses an optimization algorithm to ensure acceptable levels for the $ISPTA_{3,3}$, the Mechanical Index (MI), and the Thermal Index (TI) during use. When imaging is activated, you cannot adjust imaging parameters that affect the transmit power output until the acquisition is complete.

Quantitative Ultrasound Technology

The following features use quantitative ultrasound (QUS) technology to obtain qualitative visualizations or quantitative measurements of tissue stiffness:

- Virtual Touch UDFF (ultrasound-derived fat fraction)

QUS technology measures the attenuation coefficient (AC), backscatter coefficient (BSC), and the ultrasound-derived fat fraction index (UDFF).

Virtual Touch Point Shear Wave Elastography (pSWE)

Use Point Shear Wave Elastography (pSWE) to measure tissue for a selected region of interest (ROI) as shear velocity (Vs) or Young's Modulus of Elasticity (E).

Note: Shear waves do not propagate through fluid; do not use pSWE to measure fluids.

Note: The ultrasound system assigns a label to each acquisition. You can change the label, if necessary.

To activate Point Shear Wave Elastography (pSWE):

1. Tap **Virtual Touch** during 2D imaging and then tap **pSWE**, if necessary.
A region of interest (ROI) displays on the image.

2. To reposition the ROI, roll the trackball.

Note: You cannot adjust the size of the ROI, and the software restricts the depth of the ROI.

3. To display the measured results for elasticity on the image screen, tap **kPa**.

4. To begin acquisition, press **UPDATE**.

An audible tone indicates when pSWE measurement acquisition ends.

When acquisition is complete, the system automatically freezes, selects the active label, and displays the tissue shear velocity (Vs), elasticity (E), and the depth of the ROI in the measured results.

During system freeze, a cooling timer indicates the remaining time (in seconds) until the transducer is available for image acquisition.

Note: If the system displays the measurement values as **X.XX**, the confidence interval threshold for measurement quality was not reached. Repeat the acquisition until the system displays a numerical value instead of **X.XX**.

5. Press **SET** to confirm the acquisition.
6. To reactivate pSWE after the automatic freeze period ends:
 - a. Press **FREEZE** to unfreeze the image.
 - b. Adjust the ROI position, if necessary.
 - c. Press **UPDATE** to begin acquisition.
7. To exit, tap **Exit**.

Virtual Touch Ultrasound-Derived Fat Fraction (UDFF)

Ultrasound-Derived Fat Fraction (UDFF) measures tissue in a selected region of interest and calculates the UDFF index for assessing steatotic liver disease.

To obtain the Ultrasound-Derived Fat Fraction (UDFF):

Prerequisite: Knowledge of pSWE imaging for evaluation of tissue stiffness.

1. Tap **Virtual Touch** during 2D imaging and then tap **UDFF**, if necessary.
A region of interest (ROI) displays on the image.

2. To reposition the marker and the ROI, roll the trackball.

Note: You cannot adjust the size of the ROI.

3. Press **UPDATE** to begin acquisition.

An audible tone indicates when measurement acquisition ends.

When acquisition is complete, the ultrasound system freezes the image and displays the following values for the active measurement label in the measured results:

- Percentage of Ultrasound-Derived Fat Fraction (UDFF)
- Depth of the ROI

Note: If the system displays the measurement values as **XX**, the confidence interval threshold for measurement quality was not reached. Repeat the acquisition until the system displays a numerical value instead of **XX**.

4. Press **SET** to confirm the acquisition.
5. To obtain another UDFF measurement after the automatic freeze period ends, unfreeze the image and repeat steps 2 through 4.
6. To exit, tap **Exit**.

Measurement Labels

See also: For information about the measurement function, refer to Chapter C1 in this manual.

Point Shear Wave Elastography (pSWE) and Ultrasound-Derived Fat Fraction (UDFF) provide a set of shear velocity and elasticity measurement labels for abdomen exams:

- Up to ten sites for generic locations
- Up to ten masses
- Up to nine liver segments based on the Couinaud segmentation system

After each measurement, the system displays a counter below the image to indicate the number of valid measurements and the total number of all measurements assigned to each liver assessment measurement label. A valid measurement is any measurement that displays a numerical value in the measured results.

Site 1: 8 Valid / 10 Total

Example of measurement counter.

Note: You can make as many measurements as required, however, the ultrasound system can only save up to 20 measurements. When you make the twenty-first measurement, the ultrasound system deletes the earliest measurement and saves the latest measurement; the first measurement is replaced with the twenty-first measurement.

To change the assigned measurement label:

1. Tap **Site**, **Liver Segment**, or **Mass**.
2. Tap the required label from the list.
3. To save the updated label to the patient report, press **SET**.

Virtual Touch Patient Report and Left Panel

When you have made Point Shear Wave Elastography (pSWE) or Ultrasound-Derived Fat Fraction (UDFF) measurements, the abdomen patient report lists measurements for each liver assessment measurement label. The abdomen report also includes an assessment section with specific selections for the liver.

During imaging, the left panel displays the completed measurements and statistics for each measurement label and includes overall statistics for the measured anatomy.

To view results in the left panel:

- Click **pSWE** to display pSWE results.
- Click **UDFF** to display UDFF results.

Virtual Touch Measurement	Description
Vs (m/s)	(Available only for Point Shear Wave Elastography measurements) Shear velocity
E (kPa)	(Available only for Point Shear Wave Elastography measurements) Elasticity
UDFF (%)	(Available only for Ultrasound-Derived Fat Fraction measurements) Ultrasound-Derived Fat Fraction index for the region of interest
Depth (cm)	Depth of the region of interest
Median	Median calculation for all measurements associated with the assigned label
Mean	Mean calculation for all measurements associated with the assigned label
Std Dev	Standard deviation calculation for all measurements associated with the assigned label
IQR	Interquartile range (the distance between the 75th percentile and the 25th percentile) for all measurements associated with the assigned label
IQR/Median	Ratio of the interquartile range to the median calculation for all measurements associated with the assigned label
Overall Statistics	---
Median	Median calculation for all measurements included in the report
Mean	Mean calculation for all measurements included in the report
Std Dev	Standard deviation calculation for all measurements included in the report
IQR	Interquartile range calculation for all measurements included in the report
IQR/Median	Ratio of the interquartile range to the median calculation for all measurements included in the report

To delete a measurement:

- Click **Delete**.



Delete icon.

Scanning Techniques to Mitigate Technical Limitations

- To obtain reproducible results in the liver with Virtual Touch applications, take measurements in the right lobe of the liver using an intercostal scanning approach. Advise the patient to breathe normally and to momentarily stop breathing and then take measurements. Avoid taking measurements during deep inspiration and breath hold as the literature cites that this action increases central venous pressure. Elevated central venous pressure can artificially elevate shear wave velocity and elasticity measurements.


Siemens Healthineers recommends maintaining a measurement region of interest (ROI) depth between 3 cm and 6 cm to minimize the effects of depth on shear wave velocity and elasticity measurements.

- To obtain reliable and repeatable results with Virtual Touch applications, apply the appropriate amount of scanning pressure.

Siemens Healthineers recommends using minimal to mild compression for optimal results. Applying excessive pressure with the transducer on the body may artificially elevate shear wave velocities and result in unreliable shear wave velocity and elasticity.

- In general, use a transducer angle of 90 degrees to the skin surface for the highest measurement accuracy.


An out-of-plane transducer angle less than 50 degrees to the skin surface can result in artificially low shear wave velocities and elasticity values due to loss of transducer contact with the flat scanning surface.

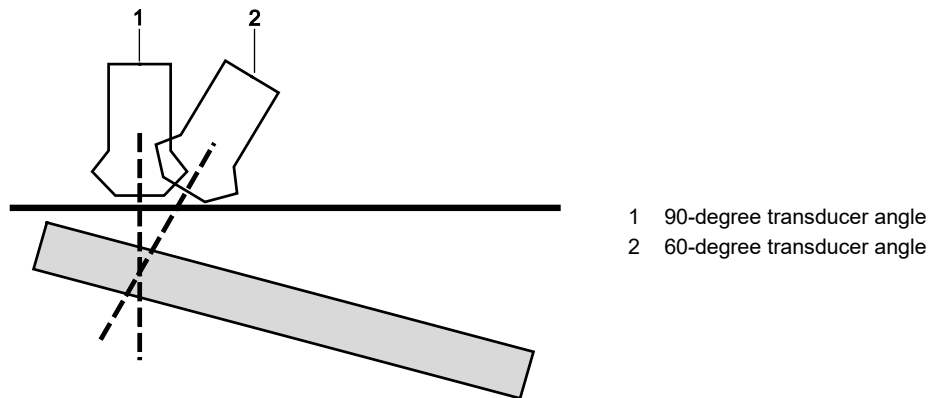

Transducer angle on a flat surface.

- Maintain adequate contact and be aware of the effects of the transducer angle relative to the contact surface.

When the transducer surface is in continuous contact with the skin surface and follows the contour of the curved surface, shear wave velocity and elasticity measurements are accurate and stable.

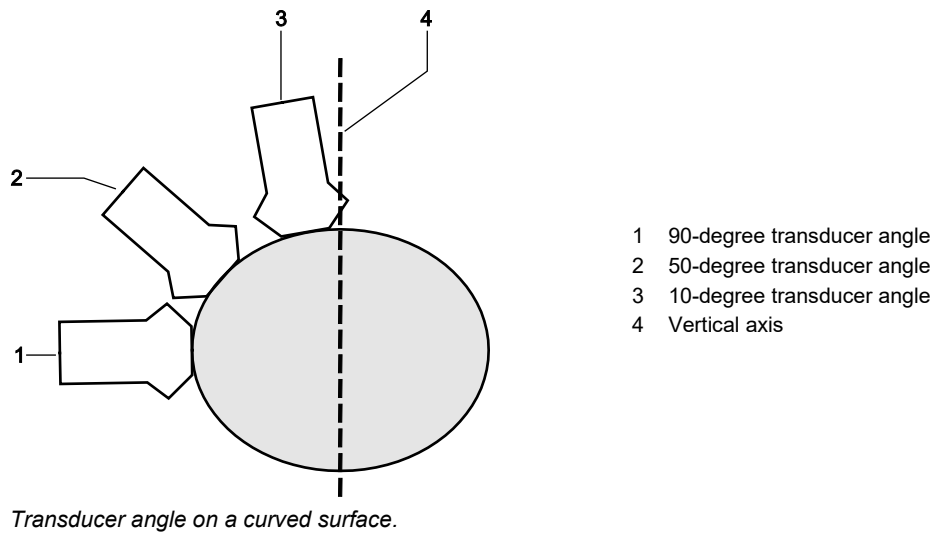

Shear Wave Velocity and Elasticity Measurements in Liver Assessment

Virtual Touch applications can provide a reliable measure of shear wave velocity and elasticity in the liver when proper technique is used.

The following factors can result in abnormally high values in the liver:

- Excessive transducer pressure on the liver, such as in subcostal scanning where transducer pressure is used to displace intervening bowel gas
- Elevated central venous pressure due to heart failure or valsalva maneuvers (including deep inspiration and breath hold)
- Significant portal hypertension
- Poor contact with the skin surface or oblique transducer angles relative to the liver
- Elevated aminotransferase levels
- Food intake less than 2 to 3 hours prior to examination
- Liver location (liver segment measured)
- Liver fibrosis

The following techniques should be used in liver assessment to ensure measurements using Virtual Touch applications in liver assessment are reliable:

- Measurements are most reliable at least 1 cm to 2 cm deep to the liver capsule using the intercostal scanning approach along a mid-axial line in liver segments 5 and/or 8.
- The patient should breathe normally with measurements taken after the patient is asked to momentarily suspend respiration.
- Multiple measurements in the same location should be taken to ensure measurement repeatability with the median shear wave velocity and elasticity value representing the most reliable single measure when there are measurement outliers.
- The Virtual Touch applications measurement region of interest (ROI) should be maintained between 3 cm and 6 cm depth whenever possible to minimize the effect of depth on measurements. Values can be higher at depths less than 3 cm and lower at depths greater than 6 cm.
- An IQR/median shear wave velocity ratio less than 0.15 or an IQR/median shear wave elasticity ratio less than 0.3 is indicative of reliable liver assessment with Virtual Touch applications. A higher ratio indicates significant variability in the measurements, decreasing the reliability of measurement results.

A7 Contrast Agent Imaging

Overview.....	2
Using Contrast Agent Imaging	3
Techniques for Contrast Agent Destruction	5
Configuring the Flash Sequence	5
Basic Measurements	6

Overview

 **WARNING:** Carefully follow the contrast agent manufacturer's instructions for use, including indications and contraindications, when administering ultrasound contrast agents.

 **WARNING:** Cardiac rhythm disturbances during perfusion studies using gas ultrasound contrast agents have been observed in the diagnostic range of the Mechanical Index (MI) values. For details, refer to the specific package insert for the contrast agent being used.

Contrast agent imaging is an enhanced form of ultrasound imaging with a wide-band, imaging technology designed for use with contrast agents.

Contrast agent imaging includes the following acquisition techniques:

- Contrast Harmonic Imaging (CHI): Amplifies the second harmonic detected while suppressing the signals from the tissue or agent within the fundamental frequency bandwidth.
- Contrast Pulse Sequence (CPS): Enhances the specificity and sensitivity of the non-linear signals from the tissue or agent within the fundamental frequency bandwidth.

Contrast agent imaging provides preconfigured settings for filters, overall gain, and other parameters coordinated with Siemens Healthineers patented data sequencing techniques for optimal viewing of contrast agent.

During contrast agent imaging, the system displays the mechanical index (MI) and the thermal indices (TI) in the imaging parameters.

To activate contrast agent imaging:

1. Tap **Contrast**.
The system activates a single display format for the cardiac exam and a dual display format for other supported exams.
2. To exit, tap **Exit**.

Using Contrast Agent Imaging

Note: Contrast agent imaging supports pulsed wave Doppler and continuous wave Doppler for the cardiac exam.

Prerequisite: Activate contrast agent imaging.

To	Do This
Enhance areas of low contrast agent concentration	(Not available for the cardiac exam or during the burst process) ▪ Tap Peak Hold to enable maximum intensity projection.
Adjust transmit power	(Not available during the burst process) ▪ Rotate [Tx Power].
Adjust the frame rate	(Not available for the cardiac exam or during the burst process) ▪ Rotate [Frame Rate].
Increase or decrease brightness of the displayed image	▪ During the cardiac exam, rotate 2D. ▪ During other supported exams, rotate [Contrast Gain].
Identify the ultrasound contrast agent	(Not available for the cardiac exam or during the burst process) ▪ Tap Contrast Agent and then select a setting. – Low MI to optimize the settings for low-range MI contrast agents. – Mid MI to optimize the settings for mid-range MI contrast agents.
Start the timer	1. Tap Timer. The system displays the elapsed time on the image screen. 2. To stop and remove the timer from the image screen, tap Timer again.
Display the recorded time during CINE	▪ Press FREEZE.

To	Do This
Emphasize the contrast agent, the tissue, or both	(Not available for the cardiac exam) 1. Tap the required setting for the left image of a dual display format. – Contrast Only to emphasize the contrast agent. – 2D Only to emphasize the tissue. – 2D + Contrast to combine of the contrast agent and tissue. 2. To display the left image in full-screen format, tap Live Dual. 3. To redisplay the image in dual format, tap Live Dual again.
Change the direction for storing a clip	(Available only for the cardiac exam) 1. Tap Clip Settings. 2. Tap Capture Type and then select a setting. – Prospective acquires a clip of the newly acquired images. – Retrospective acquires a clip of the previously stored images.
Change the method for determining the duration of the clip	1. Tap Clip Settings and then select a setting. – Time defines the length of the clip to capture in seconds. – Beat defines the length of the clip to capture in heart cycles. – R-Wave defines the capture point based on a delay beyond the R-wave for cardiac exams. 2. Tap the arrow to increase or decrease the clip duration.

Techniques for Contrast Agent Destruction

You can use a high mechanical index technique to destroy the contrast agent microbubbles and then observe reperfusion of the bubbles into the anatomy.

Contrast agent imaging includes the following techniques for contrast agent destruction:

- Burst is a high mechanical index technique for rapid destruction of contrast agent microbubbles for observation of the time required for reperfusion of the contrast agent.
- Flash sequencing is a high mechanical index non-imaging mode technique for rapid destruction of contrast agent and observation of the reperfusion of the contrast agent into the anatomy.

Note: Flash sequencing is not available for the cardiac exam.

To use the Burst destruction technique:

1. To specify the duration for destruction of the contrast agent, rotate **[Burst Time]**.
2. Tap **Burst**.
3. To stop the burst process before the specified burst time ends, tap **Burst** again.

To use the Flash Sequence technique:

1. Configure the Flash Sequence, if necessary.
2. Tap **Flash Sequence**.
The system activates the continuous burst of the contrast agent destruction and reperfusion sequences.
3. To cancel the flash sequence, tap **Flash Sequence** again.

Configuring the Flash Sequence

You can configure the flash sequence using selections on the touch screen.

To configure the Flash Sequence:

1. Tap **Configure Flash Sequence**.
2. Select the first step of the sequence.
 - a. Tap **Burst Time**.
 - b. Tap the duration, in seconds, for destruction of the contrast agent.
3. Select the second step of the sequence; tap an option.
 - **Peak Hold** starts maximum intensity projection when Burst is complete.
 - **Live Imaging** resumes real-time imaging when Burst is complete.
4. Tap **Imaging Duration** and then select the length of acquisition in seconds.
5. Select the third step of the sequence; tap an option.
 - **Stay Live** continues real-time imaging when acquisition is complete.
 - **Auto Freeze** freezes the image when acquisition is complete.
6. To capture the flash sequence as a clip, tap **Flash Sequence Clip Capture**.
7. When configuration is complete, tap **Exit**.

Basic Measurements

You can perform all basic measurements during contrast agent imaging on either image in a dual screen layout or simultaneously on both images using the shadow function.

See also: For information on measurements and calculations, refer to Chapter C1 in the Clinical Applications Guide.

To perform a measurement using the shadow function:

1. Activate the measurement function.
2. Tap the required measurement tool.
3. Tap **Shadow**, if necessary.
4. Measure the area of interest on either image.
The ultrasound system duplicates and displays the measurement on the other image.
5. To exit the shadow function, tap **Shadow**.

B1 Physiologic Function

Physiologic Function	2
ECG Function	2
Connecting ECG or Auxiliary Cables	3
Connecting Leads to the Patient	4
Using the Auxiliary Cable	5
Activating the ECG or AUX ECG Function	7
Adjusting the Traces	7
Working with the Trigger Function	8
Displaying Auxiliary Signals	8
ECG Troubleshooting Guide	9

Physiologic Function

The physiologic function includes the ECG function and physiologic inputs. The system can display auxiliary physiologic traces from medically approved, properly isolated devices. You can simultaneously display an ECG trace and a respiratory trace.

ECG Function

 **WARNING:** Do not use the ECG option in conjunction with electrosurgery or diathermy equipment.

 **WARNING:** Use ECG for timing purposes only. It is not intended for diagnostic usage or patient monitoring.

Electromagnetic Compatibility (EMC) Note: Operating the ultrasound imaging system in close proximity to sources of strong electromagnetic fields, such as radio transmitter stations or similar installations, as well as portable and mobile RF (radiated frequency) may lead to temporary degradation or interference visible on the monitor screen. However, the device has been designed and tested to withstand such interference and will not be permanently damaged. Wait until the interference disappears or restart the ultrasound system.


Note: If the physio module detects abnormal input due to electrical interference, the module is designed to shut down to prevent the incorrect display of the physio signals on the image screen. If the physio module shuts down, you must restart the ultrasound system to reset the physio module.

Note: The ECG Inputs are defibrillation proof. However, in the event of defibrillation during the ECG function, the ECG inputs may become saturated (overloaded). An ECG pattern may not be visible for up to 30 seconds. After 30 seconds, the ECG trace displays on the screen.

Connecting ECG or Auxiliary Cables


⚠ WARNING: To reduce risk of electric shock and burns, use only the cable and patient leads supplied with the physiologic function. Use of other ECG cables could defeat the current-limiting and electrical safety features of the physiologic function.

⚠ Caution: To avoid possible damage to the physio cables, do not pull on the cables when disconnecting the ECG connector. Always grasp the connector housing.

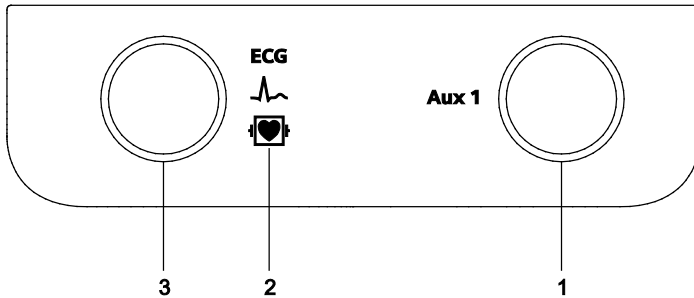

Example of physio panel on the left side of the ultrasound system.

- 1 Auxiliary port
- 2 ECG port label
- 3 ECG port

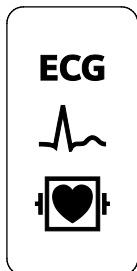

Example of ECG port label.

To connect the ECG cable:

- Connect the 6-pin, 3-lead ECG cable to the ECG port on the ultrasound system.

To connect the Auxiliary cable:

- Connect the 7-pin, auxiliary cable to the auxiliary port on the ultrasound system.

Connecting Leads to the Patient

In keeping with existing international standards, the ECG leads for systems that operate at 115V, for example North and South America, are red, white, and black. Leads for systems that operate at 230V, for example Europe, are green, red, and yellow. Use only ECG leads that comply with the requirements of your country.

The system voltage corresponds to the following ECG leads:

System Voltage	Green	Red	Yellow	White	Black
100V	Left leg	Right arm	Left arm	--	--
115V	--	Left leg (LL)	--	Right arm (RA)	Left arm (LA)
230V	Left leg (F)	Right arm (R)	Left arm (L)	--	--

To connect the leads to the patient:

1. Attach the ECG electrode patches to the patient in the locations specified above.
2. Connect each lead to the appropriate patch.

Use the configuration settings to select the lead type for tracing the voltage between the electrodes (right arm and left arm or right arm and left leg).

Using the Auxiliary Cable

⚠ WARNING: To reduce the risk of electrical shock while using the auxiliary port on the ultrasound system, ensure that the external source (peripheral equipment) is properly set up and is designed for direct patient connection.


The stress adapter cable kit contains a stress adapter cable and connector adapters used to attach various peripheral equipment to the ultrasound system.

See also: For information about connecting to peripheral equipment, refer to Chapter 3 in the Instructions for Use.

To connect an auxiliary cable to the ultrasound system and to the peripheral equipment:

1. Attach the auxiliary input connector to the auxiliary adapter.


Note: The auxiliary connector connects to the ultrasound system.

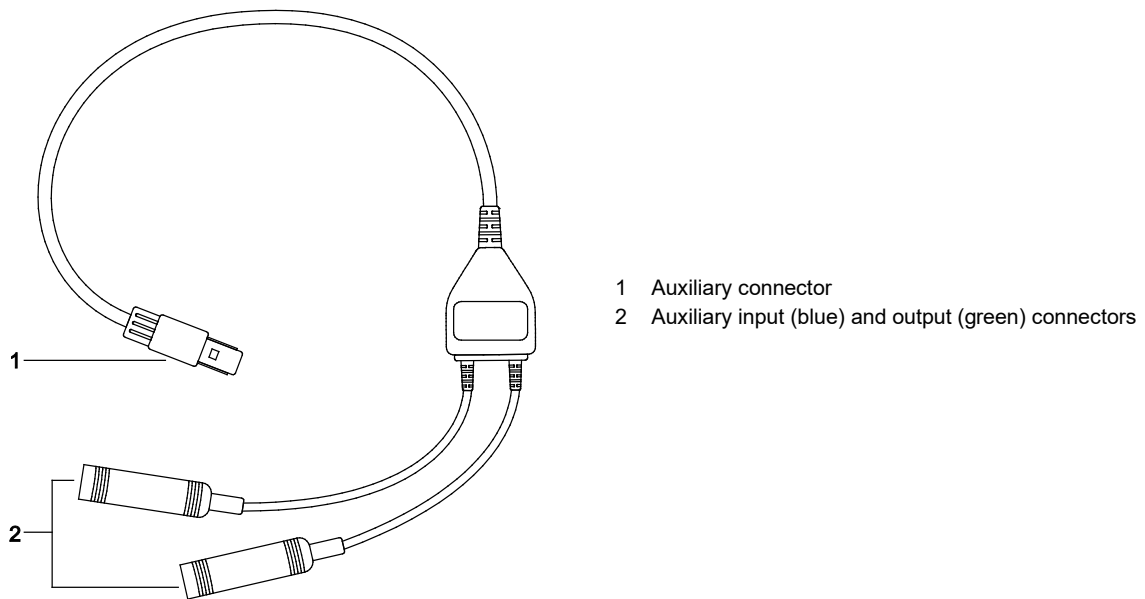


2. Attach one end of the stress adapter cable to the auxiliary adapter.

3. Attach the required connector adapter for the peripheral equipment to the other end of the cable.


Note: Refer to the manufacturer's operating instructions for the peripheral equipment to determine the appropriate connector adapter to use.

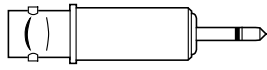

Phone adapter.

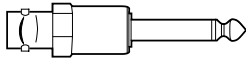

Phone plug (mono).

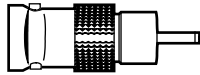

BNC to Phono plug adapter.

Activating the ECG or AUX ECG Function

You can enable the system to display an ECG trace, an auxiliary trace, or an R-wave indicator, if specified in the configuration settings.

To display an ECG trace:

Note: If baseline wander in the ECG trace becomes a problem, the cause may be respiratory movement (ventilatory excursions). Minimize this artifact by making certain that the ECG leads or cable are not needlessly draped over the patient. Secure the leads and cable off to the side of the patient.

1. Connect the ECG cable to the ECG port.
2. Connect the ECG leads to the patient.
3. Tap **Physio Settings** and then tap **ECG**.
During CINE, the system displays a yellow vertical line to indicate the relationship of the specific phase of the ECG trace and two purple markers to indicate end diastole and end systole.
4. To activate the respiratory trace during an ECG trace, tap **Resp**.
5. To exit the physio function, tap **ECG** and then tap **Resp**.

To display an auxiliary trace:

1. Connect the auxiliary cable to the auxiliary port.
2. Connect the auxiliary cable to the external ECG monitoring device.
3. Tap **Physio Settings** and then tap **Aux ECG**.
During CINE, the system displays a yellow vertical line to indicate the relationship of the specific phase of the ECG trace and two purple markers to indicate end diastole and end systole.
4. To exit the physio function, tap **Aux ECG**.

Adjusting the Traces

You can adjust the trace gain and trace position during the exam.

To adjust the trace gain and position:

1. To adjust the gain of the trace, rotate **[Gain]**.
2. To adjust the vertical position of the trace, rotate **[Position]**.
3. To invert the trace, tap **Invert**.
4. To activate automatic gain control, tap **AGC**.

Working with the Trigger Function

You can use triggers at selected points in the cardiac cycle when the system refreshes the 2D image. The trigger function is not compatible with the following modes or features:

- M-mode
- Doppler
- Compounding

To enable single or dual triggers:

- To display a single trigger, tap **Trigger**.
- To display dual triggers, tap **Trigger** and then tap **Trigger Type**.

The system detects the R-wave and displays each trigger as a white, vertical line on the ECG trace. During Cine, the system displays a yellow vertical line to indicate the relationship of the specific phase of the ECG trace to the 2D image.

To set the trigger delay:

1. Rotate **[Ch 1 Delay]**.

Note: The position of the trigger advances in increments of ten milliseconds.

Note: When you change a trigger delay, a lag period occurs before the new trigger position is activated.

2. To position the second trigger, rotate **[Ch 2 Delay]**.

Note: If you select two triggers and the triggers have identical values, then the system derives the second update time from the frame rate.

Displaying Auxiliary Signals

You can display auxiliary signals such as respiratory traces or R-wave detection on the ultrasound system.

To display auxiliary signals:

- To display auxiliary signals on the system, tap **Physio Settings**, and then tap **Aux Signal**.

ECG Troubleshooting Guide

Symptom	Possible Cause
ECG trace is not present, or cannot detect an R-wave	▪ The ECG cable is not properly attached to the system. ▪ Leads are poorly connected to the cable. ▪ Leads are poorly connected to the patient. ▪ Leads are incorrectly placed on the patient. Remove the ECG patch, clean the skin, and attach a new ECG patch to the patient. ▪ ECG gain is set too low. Rotate [Gain] to adjust the wave amplitude of the trace. Or, tap AGC to enable the automatic gain control.
Cannot position a trigger	▪ No R-wave present (see above). ▪ Too much movement by the patient.
ECG trace is noisy, ragged, or erratic	▪ ECG gain is set too high. Rotate [Gain] to adjust the wave amplitude of the trace. Or, tap AGC to enable the automatic gain control. ▪ A lead or leads may be detecting muscle movement. Inspect the ECG patch placement and reposition the lead(s) on the patient as necessary. ▪ The conductive gel pad of an ECG patch has evaporated (dried out). Replace the ECG patch with a new one (check the expiration date of the new patch).
Trigger is not visible on the trace	▪ A trigger may be superimposed on the R-wave. ▪ If dual triggers are used, the first and second trigger may be superimposed. Select Ch 1 Delay or Ch 2 Delay to reposition a trigger.

B2 Stress Echo

Overview.....	2
Working with Stress Echo.....	2
Stress Echo Protocol Status Indicators	4
Acquiring Clips in a Stress Echo Protocol.....	5
Selecting Preferred Clips for Analysis	6
Sorting and Comparing Images	8
Performing Measurements	10
Wall Motion Scoring.....	10
Stress Echo Report	10

Overview

Note: Siemens Healthineers recommends that you are familiar with stress echo examinations before attempting to acquire stress echo clips. Refer to the American Society of Echocardiography (ASE) protocol in the following publication: Schiller NB, Shah PM, Crawford M, et al. "Recommendations for Quantification of the Left Ventricle by Two-Dimensional Echocardiography." *Journal of the American Society of Echocardiography* 2(5):364, 1989.

You can capture and review cardiac clips for multiple-phase (multiple-stage) predefined and user-defined stress echo protocols. A cardiac clip displays motion from beginning systole to end systole (end diastole to end systole) or R-wave to R-wave, as indicated by the R-wave of the ECG trace.

Working with Stress Echo

Stress Echo acquires clips in views and stages using either continuous or non-continuous acquisition, depending on your configuration of each stage.


The system includes the following predefined stress echo protocols:

- **Pharmacology**
- **Pharmacology Continuous**
- **Bicycle**
- **Bicycle Continuous**
- **Treadmill**
- **Treadmill Continuous**

You can use a predefined protocol or create a user-defined protocol to configure views, stages, and settings, including prospective or retrospective acquisition, the number of clips per capture, the capture length, the capture delay, the display of the stage timer, and the display of the review screen.

When you complete acquisition for each stage, the system advances to the next task in the protocol as defined in the protocol configuration.

Use the configuration settings to create and modify stress echo protocols.

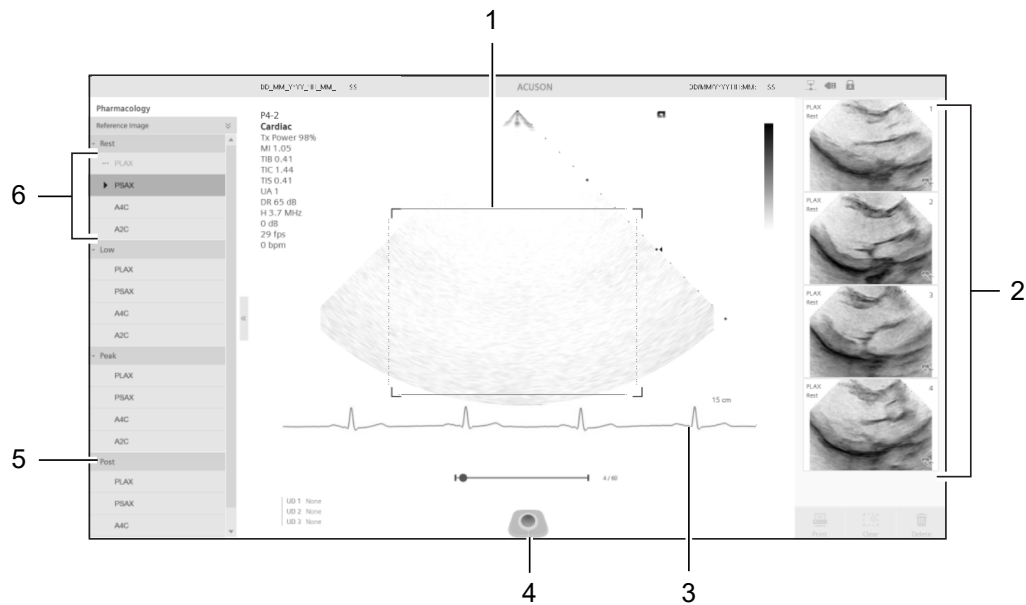

Example of the setup screen for stress echo.

- 1 Region of interest (ROI)
- 2 Image thumbnails acquired during stress echo
- 3 ECG trace
- 4 Trackball legend
- 5 Stages

The total number of acquired clips and the number of remaining captures displays next to the label for each stage during a continuous protocol.

- 6 Views

The total number of acquired clips displays next to the label for each view during a continuous protocol.

Stress Echo Protocol Status Indicators

The following icons display next to the study name in the patient browser to indicate the status of the stress echo protocol.

Icon	Description
	The stress echo protocol is activated.
	The stress echo protocol is paused.
	The stress echo protocol is completed.

To activate a stress echo protocol:

- Tap **Stress Echo** and then select the required protocol.

To	Do This
Pause a stress echo protocol	Tap Pause Protocol. An on-screen symbol indicates the protocol is paused. 
Complete a stress echo protocol	Tap End Protocol. Note: You cannot modify a stress echo protocol after the end of the protocol. Note: You can only export a completed stress echo protocol.
Restart a paused stress echo protocol from a closed study	- Access the patient browser and then select the required study. - Tap Stress Echo.
Select a different stress echo protocol	Complete the stress echo protocol and then select a new study. Note: Only one stress echo protocol is available for each study.

Acquiring Clips in a Stress Echo Protocol

Use the stress echo configuration settings to enable the timer.

To begin acquiring clips:

1. Adjust the position or size of the capture region of interest, if necessary.
Note: The adjusted depth, focus, position, and size of the region of interest in the first stage are applied to subsequent acquisitions of the same view.
2. Press **CLIP STORE**.
The system begins acquiring and labeling clips and starts the stage timer, depending on the selected protocol. If you pause a protocol, the stage timer continues to count the elapsed time.

When the capture is completed, the background of the view label changes to three dots.

During a non-continuous protocol, the system stops acquiring clips depending on the setting for capture. During a continuous protocol, the system acquires clips before you press **CLIP STORE** again.


To	Do This
Display the ECG trace	▪ Tap Physio Settings .
Activate another view or stage in a workflow protocol	▪ Tap the required view label. Or, rotate [Stage] to select the required stage label.
Start or stop the manual timer	▪ Tap Timer .
Acquire clips with color activated	1. Press C . 2. Press CLIP STORE .
Reacquire clips for the current view	▪ Select the label for the required view during acquisition.

Selecting Preferred Clips for Analysis

You can identify clips for review, comparison, and analysis. The system activates the selection task according to the protocol configuration, for example, at the end of each view, at the end of each stage, or at the end of all stages.

The selection screen displays the acquired clips for views and stages. All images and clips acquired during the exam are numbered in the sequence of acquisition.

Acquired clips not selected are removed at the end of the exam. If you do not complete the selection task before you end the stress echo exam, the system saves all acquired clips.

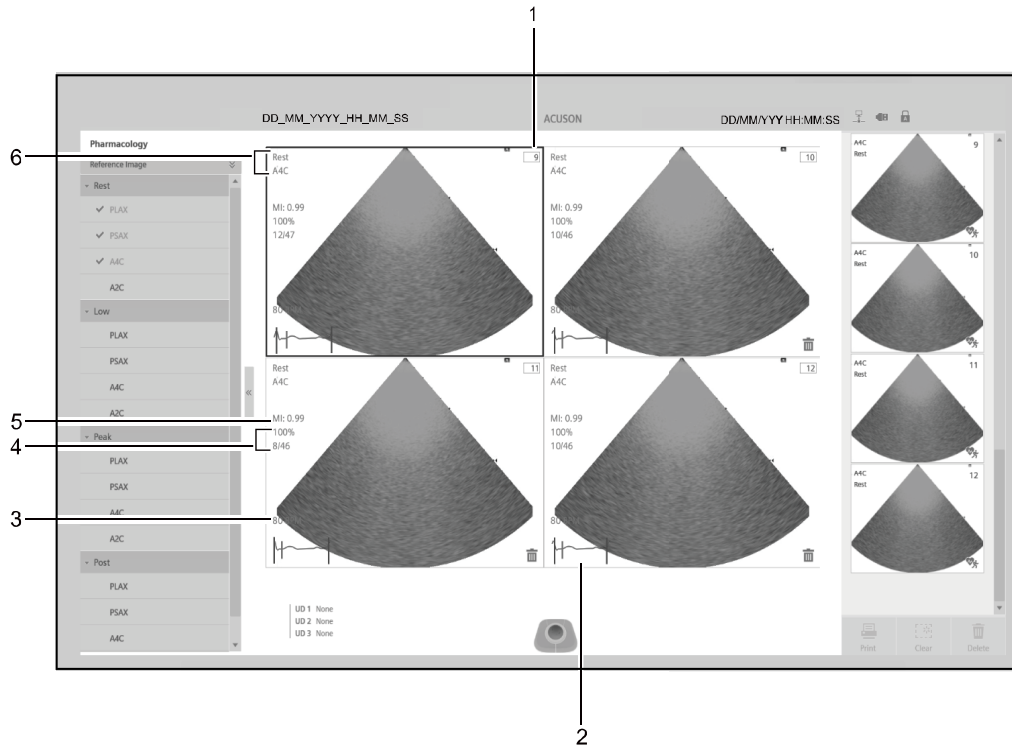

Example of the selection screen during stress echo.

- 1 Sequence of clips
- 2 ECG trace
- 3 Heart rate
- 4 Playback speed and number of frames indicated by the current frame number and the total frame number
- 5 Mechanical index
- 6 Captured stage and view

To select preferred clips after capturing a view or stage:

1. Tap **Accept All**.
The system displays a blue border around the acquired clips and saves the acquired clips. The system closes the selection window and activates the next view or stage in the protocol.
2. To postpone the selection task, tap **Defer**.
The system closes the selection window and activates the next view or stage in the protocol.
3. To reject all currently acquired clips and reacquire clips for the current view, tap **Reject All**.
The system closes the selection window and reactivates acquisition for the current view.

To select preferred clips during stress echo review:

1. Click the preferred clips.
The system displays a blue border around the selected clips and saves the selected clips.
2. To deselect a preferred clip, click the clip.
The system removes the blue border from the clip and displays the delete icon to indicate the clip will be removed when you complete the protocol.



Sorting and Comparing Images

During review, the system displays only the clips selected as preferred clips.

To activate stress echo review:

- Tap **Review Mode**.

To	Do This
Cycle through a set of selected clips	▪ The system displays the reverse sequence icon.  ▪ Click the forward button to view clips in a forward sequence. ▪ Click the backward button to view clips in a reverse sequence.
Sort the acquired images by view	▪ Tap Sort in the left panel and then tap By View. Stress review displays the first view for all stages, followed by each subsequent view for all stages.
Sort the acquired images by stage	▪ Tap Sort in the left panel and then tap By Stage. Stress review displays all the views for the first stage, followed by the views for each subsequent stage.
Compare views for two stages	1. Tap Compare. 2. Select the stages to compare from the Reference and Compare lists.

To activate review:

- Activate stress echo and then press **REVIEW**.

To	Do This
Display the first page of clips acquired for each view	▪ Tap Go to Next Group .
Display consecutive pages of clips	▪ Tap Page .
Match the playback synchronization of clips	▪ Tap Sync .
Select a segment of clips for playback	▪ Tap Mode and then select a setting: – Full R-R displays the full heart rate cycle of systolic and diastolic portions. – Adjusted displays the user-adjusted segment of clips. – Diastole only displays only the diastolic portion of the heart cycle. – Systole only displays only on the systolic portion of the heart cycle.
Display an end-diastolic frame or end-systolic frame	(Available only in single display format) ▪ Tap ED Frame or ES Frame .
Display an end-diastolic frame and end-systolic frame simultaneously	▪ Tap ED/ES Split .
Display or hide information, for example, the heart rate	▪ Tap Text Overlay .
Return to stress echo acquisition mode	▪ Press REVIEW .
Display the review screen without exiting a stress echo protocol	▪ Tap General Review .
Change a view label	Note: The function is only available when the selection mode is activated. ▪ Tap Rename View .

Performing Measurements

You can create a trace to measure the volume of the left ventricle (LV) for stress echo.

Note: Measure and then label is not available during a stress echo protocol.

To activate measurements:

1. Select the required clip in a single display format or ES/ED split format in the stress echo review screen.
2. Activate the measurement function.
3. Tap the label on the measurement menu.
The system activates the Simpson BP measurement labels.
4. Position and anchor all the measurement markers required for the measurement.

See also: For information about the measurement function, refer to Chapter C1 in this manual.

Wall Motion Scoring

The system provides the function for user-assigned wall motion scores. You can assign wall motion scores (WMS) to specific segments within each view (representative loop). You can also assign a normal wall motion score to the currently selected view or to all displayed views.

To assign wall motion scores:

1. Tap **WMS Scoring**.
2. To assign individual scores to each wall motion segment, select a number.
The system displays the core levels or options for scoring.

To customize the display of wall motion scoring:

1. Tap **Configuration** and then tap **Measurement**.
2. Click **General**.
3. Select a **Report View Type** from the **Wall Motion Scoring** selections, if necessary.
 - **Bull's eye**
 - **View Plane**
 - **Both**

Stress Echo Report

Note: The measurement report for stress echo is displayed only during a stress echo protocol.

To display the report page for stress echo during a stress echo protocol:

- Press **REPORT**.

C1 Measurements and Calculations

Measurement Tools and Procedures	2
Activating the Measurement Function	3
Measurement Methods.....	3
Measurement Markers.....	4
Measurement Labels	5
Measurement Results	5
Label then Measure or Measure then Label	6
Editing Measurements.....	6
General 2D-Mode Measurement Methods	7
General M-Mode Measurement Methods	13
General Doppler Measurement Methods	14
Follicle Measurements	18
Performing Pediatric Hip Measurements.....	19
Calculated Results of the Pediatric Hip Tools.....	21
Measuring Carotid Intima-Media Thickness with the Auto IMT Feature	22
Calculated Results of Auto IMT Tools	23
Patient Reports	24
Modifying Reports.....	25
Tracing Areas of Interest	26
Printing and Exporting Reports.....	26
DICOM Structured Reporting	26
Pediatric Hip Report	27

Measurement Tools and Procedures

The measurement function provides general measurements and calculations for all exam types and the following images:

- Real-time or frozen 2D images
- Frozen M-mode sweeps and Doppler spectrums
- Stored images

Note: You can make measurements during imaging or during review.

The system also provides exam-specific measurements and calculations for some exam types. Calculations require one or more specific measurements. When you complete all the measurements required by a calculation, the system displays the result.

Use the configuration settings to customize the measurement function for each exam type.

See also: For information about customizing the measurement function, refer to Chapter 3 in the System Reference.

Activating the Measurement Function

Use the configuration settings to configure the **FREEZE** key to activate the measurement function.

To	Do This
Activate the measurement function during imaging	- For 2D-mode, press CALIPER. Or, press FREEZE, if specified in the configuration settings. - For M-mode and Doppler, freeze the sweep or spectrum and then press CALIPER.
Activate the measurement function on stored images or clips	Note: Images must be stored in ultrasound format, not screen capture format. - Press REVIEW. - Select the required image or clip. - Press CALIPER. - Press the control assigned to the image store function (for example, PRINT STORE 1) to save the changes. - To store the measured results, press REPORT and then press the control assigned to the image store function (for example, PRINT STORE 1). The system saves the measured results and a screenshot of the report.
Exit the measurement function	Unfreeze the image. Or, press CALIPER.

Measurement Methods

When the measurement function is active, the measurement methods available for the active exam type and operating mode are displayed on the touch screen.

Use the configuration settings to customize the default measurement method for each mode and the default method for each type of calculated measurement, for example, volume.

To	Do This
Select a measurement method	Note: If you label first and then measure, the system selects the measurement method. Tap a measurement method.
Access additional measurement tools	Tap More. 

Measurement Markers

Measurement markers indicate the beginning and ending points of each measurement. A number or object identifies each set of measurement markers. When a measurement method requires multiple measurements, for example, volume, the system identifies all the required markers with the same number or object.

You can display a maximum of twenty sets of measurement markers on an image at one time. When twenty marker sets are displayed, you must delete one or more marker sets before you can make another measurement.

Use the configuration settings to select an identifier for each set of measurement markers and the initial position of the measurement marker when you activate the measurement function for each exam type.



Example of measurement markers

Note: To prevent loss of a measurement in progress, always anchor the final marker before you display another screen or exit the measurement function.

To	Do This
Position a measurement marker	Roll the trackball. Measurement markers ready to be placed are green.
Anchor a measurement marker	Press SET. Or, if indicated in a procedure in this manual, press SET twice. Anchored measurement markers are yellow.

Measurement Labels

When you activate the measurement function, the system displays measurement labels for the active exam type and operating mode on the measurement menu. You must label a measurement to include the results in a report.

When you select a label, the system activates the default measurement method for the label.

Use the configuration settings to create or delete user-defined labels, to specify the labels to display on the measurement menu and in reports, and to arrange the sequence of labels on the menu.

To	Do This
Select a measurement label	1. Tap a prefix on the measurement menu, if necessary. – Right indicates the right side of the anatomy. – Left indicates the left side of the anatomy. – Prox indicates a proximal view of the anatomy. – Dist indicates a distal view of the anatomy. – Mid indicates a middle view of the anatomy. – None removes any modifier or prefix from the selected label. 2. Tap the measurement label on the measurement menu.
Display a different page of measurement labels	▪ Swipe the touch screen to view more measurement labels.
Select a label from a different exam type on the measurement menu	1. Tap a different exam type from the exam type list above the measurement menu. 2. Tap an available measurement label on the measurement menu. The system creates a report for the exam type you used for measurements.

Measurement Results

As you make measurements, the system displays the results of incomplete measurements in green and the results of completed measurements in white.

Use the configuration settings to define the default position for the measurements results, to change the background for the results, and to specify the results to display on the measurement screen and in reports.

To	Do This
Change the location of measurement results during the measurement function	1. Position the cursor on the measurement results. 2. Press and hold SET and then roll the trackball to the required location.
Change the font size of the measured results	▪ Rotate [Font Size].

Label then Measure or Measure then Label

To	Do This
Label and then measure	1. Activate the measurement function. 2. Tap the label on the measurement menu. The system activates the default measurement method for the label. 3. To use a different measurement method, tap a different measurement method soft key. 4. Position and anchor all the measurement markers required for the measurement.
Measure and then label	Note: If you measure and then select a label incompatible with the measurement method, the system retains the first measurement as an unlabeled measurement and starts a new, labeled measurement using the default method for the label you selected. 1. Activate the measurement function. 2. Tap a measurement method. 3. Position and anchor all the measurement markers required for the calculation. 4. Tap a label in the measurement menu.

Editing Measurements

To	Do This
Delete a measurement	Note: Measurement results are not removed from the report when you delete a labeled measurement from the image screen. 1. Press [Delete] to delete the most recent measurement. 2. Press [Delete] again to continue deleting measurements.
Delete all measurements from the image	Note: Measurement results are not removed from the report when you delete all measurements. ▪ Press [Delete All].
Edit a completed measurement	1. Press UPDATE to select and reposition the ending marker. 2. Press [Select Next] repeatedly to select and reposition markers in the reverse sequence.
Edit a completed trace	1. Press UPDATE to select the trace. 2. Rotate [Undo/Redo] to remove each incorrect segment. 3. Reposition and anchor the marker.
Edit an incomplete trace measurement	▪ Rotate [Undo/Redo] to remove each segment.

See also: Modifying Reports, page C1-25

General 2D-Mode Measurement Methods

Use the configuration settings to customize the default measurement method for each operating mode and the methods used for calculations.

You can make volume measurements on single or dual 2D images. If a volume calculation requires you to make measurements on multiple planes, you can use dual images, or you can change the image plane and measure single images.

Prerequisite: The measurement function is active during 2D-mode.

To	Do This
Measure distance	1. Tap Distance.  2. Position and anchor the beginning and ending markers.
Measure multiple distances	1. To measure three distances, tap 4 Point Dist.  2. To measure four distances, tap 5 Point Dist.  3. To measure five distances, tap 6 Point Dist.  4. Position and anchor the beginning and ending markers.
Determine an average distance using two measurements	1. Tap 2 Distance Avg.  2. Position and anchor the beginning and ending markers.
Determine an average distance using three measurements	1. Tap 3 Distance Avg.  2. Position and anchor the beginning and ending markers.
Measure the length of a curved structure	1. Tap Open Spline.  2. Position and anchor the first marker to start a trace. 3. Trace the structure. a. Position and anchor three or more additional markers. b. Press SET twice to anchor the ending marker.

To	Do This
Measure the circumference and area using an ellipse	Use the Trace method when the ratio of the axis lengths is more than 2:1. Use the Ellipse method when the two axis lengths are close to equal. 1. Tap Ellipse.  2. Position and anchor the beginning marker to indicate the first axis of the ellipse. 3. Roll the trackball to rotate the ellipse and size the first axis to fit the ellipse to the structure and then anchor the marker. 4. Adjust the length of the second axis to fit the ellipse to the structure and then anchor the marker.
Measure the circumference and area of an irregular shape	1. Tap Trace.  2. Position and anchor the first marker to start a trace. 3. Trace the structure and then anchor the ending marker. If the ending marker does not meet the beginning marker, the system completes the trace.
Measure the circumference and area of an irregular shape using a system-generated trace	1. Tap eSieCalcs.  2. Position and anchor the first marker to start a trace. The system displays a region of interest with a center marker. 3. Roll the trackball to align the marker centered within the rectangle to the structure of interest and then press SET twice. The system traces the structure and then activates the spline measurement tool.
Measure the length of a curved structure and area	1. Tap Spline.  2. Position and anchor the first marker to start a trace. 3. Trace the structure. a. Position and anchor three or more additional markers. b. Press SET twice to anchor the ending marker.
Measure the length of an irregular line	1. Tap Open Trace.  2. Position and anchor the first marker to start a trace. 3. Trace the structure and then anchor the ending marker.

To	Do This
Measure two angles	1. Tap Cross Angle.  2. Position and anchor the beginning and ending markers for the reference line. 3. Position and anchor the beginning and ending markers for the intersecting line to create the angle.
Measure a single angle	1. Tap Single Angle.  2. Position and anchor the beginning and ending markers for the reference line. 3. Position the ending marker.
Measure a cardiac axis angle	1. Tap Cardiac Axis.  2. Position and anchor the beginning and ending markers for the reference line. 3. Position and anchor the beginning and ending markers for the intersecting line to create the angle.
Determine a volume using one distance	1. Tap Vol 1 Dist.  2. Position and anchor the beginning and ending markers for the length.
Determine a volume using two distances	1. Tap Vol 2 Dist.  2. Position and anchor the beginning and ending markers for the length. 3. Position and anchor the beginning and ending markers for depth.
Determine a volume using three distances	1. Tap Vol 3 Dist.  2. Position and anchor the beginning and ending markers for the length. 3. Position and anchor the beginning and ending markers for the depth. 4. Position and anchor the beginning and ending markers for the width. 5. If using single 2D images, unfreeze the image, acquire and freeze the second image, and then reactivate the measurement function. The system retains the measurements from the first image when you unfreeze the image. 6. Position and anchor the beginning and ending markers for the width on the second image. The system uses the length and depth from the first image and the width from the second image to determine the volume.

To	Do This
Determine a volume using an ellipse	Note: Perform an area measurement and identify the axis common to both planes. 1. Tap 1 Ellipse.  2. Position and anchor the first marker on the axis common to both planes. The dotted line of the ellipse represents the axis common to both planes. 3. Rotate and size the axis and then anchor the ending marker. 4. Adjust the shape of the ellipse along the second axis, and anchor the ending marker.
Determine a volume using an ellipse and a distance	Make the ellipse measurement in one plane and then make the distance measurement for length in a different plane. 1. Tap Vol 1E + 1D.  2. Position and anchor the first marker on the axis common to both planes. The dotted line of the ellipse represents the axis common to both planes. 3. Rotate and size the axis and then anchor the ending marker. 4. Adjust the shape of the ellipse along the second axis, and anchor the ending marker. 5. If using single 2D images, unfreeze the image, acquire the second image in another plane, freeze the image and then reactivate the measurement function. The system retains the measurements from the first image even when you unfreeze the image. 6. Position and anchor the beginning and ending markers for the distance measurement.
Determine a volume using the method of disks trace	1. Tap Disk Trace.  2. Position and anchor the first marker to start a trace. 3. Trace the structure. a. Position and anchor three or more additional markers. b. Press SET twice to anchor the ending marker. If your ending marker does not meet the beginning marker, the system completes the trace and positions the axis line halfway between the two markers.

To	Do This
Determine a volume using the method of disks	1. Tap Disk.  2. Position and anchor the first marker to start a trace. 3. Trace the structure. a. Position and anchor three or more additional markers. b. Press SET twice to anchor the ending marker.
Determine a percentage stenosis using a distance	Note: Requires two diameters of the same vessel. 1. Tap D-%Stenosis.  2. Position and anchor the beginning and ending markers to indicate the first diameter of the vessel. 3. Position and anchor the beginning and ending markers to indicate the second diameter of the vessel.
Determine a percentage stenosis using an area	Note: Requires two cross-sectional measurements of the same vessel. 1. Tap A-%Stenosis.  2. Position and anchor the beginning marker to indicate the first axis of the ellipse. 3. Roll the trackball to rotate the ellipse and size the first axis to fit the ellipse to the vessel and then anchor the marker. 4. Adjust the length of the second axis to fit the ellipse to the vessel and then anchor the marker. 5. Repeat steps 2-4 for the second measurement.
Determine a percentage stenosis using an ellipse trace	1. Tap ET-%Stenosis.  2. Position and anchor the beginning and ending markers to indicate the first ellipse of the vessel. 3. Position and anchor the beginning and ending markers to indicate the second ellipse of the vessel.

To	Do This
Determine a percentage stenosis using auto stenosis	1. Tap Auto Stenosis.  2. Trace the structure. a. Position and anchor three or more additional markers. b. Press SET twice to anchor the ending marker.
Determine a distance using an O marker	1. Tap O-Marker.  2. Rotate SELECT to determine the O-marker size.
Determine a distance ratio	1. Tap Dist Ratio.  2. Position and anchor markers for two distance measurements.
Determine a volume ratio	1. Tap Vol Ratio.  2. Position and anchor markers for volume measurements.
Determine an area ratio	1. Tap Area Ratio.  2. Position and anchor markers for area measurements.
Determine a ratio of two dimension measurements	1. Tap d:D Ratio.  2. Position and anchor markers to measure the two dimensions and the ratio.

General M-Mode Measurement Methods

Use the configuration settings to customize the default measurement method for M-mode and the method used to determine M-mode ratios.

The horizontal axis of the M-mode sweep represents time and the vertical axis represents distance.

Prerequisite: The measurement function is active during M-mode.

To	Do This
Measure distance	1. Tap Distance .  2. Position and anchor the first marker. 3. Position and anchor the second marker.
Determine the change in distance over time	1. Tap Slope .  2. Position and anchor the first marker. 3. Position and anchor the second marker. A dotted line indicates the slope.
Determine heart rate	1. Tap HR .  2. Position and anchor the first marker at the beginning of the cardiac cycle. 3. Position and anchor the second marker at the end of the cardiac cycle.
Measure a time interval	1. Tap Time .  2. Position and anchor the beginning and ending markers.
Determine the ratio of two distances	1. Tap Dist Ratio .  2. Position and anchor markers for two distance measurements.
Determine the ratio of two time measurements	1. Tap Time Ratio .  2. Position and anchor markers for two time measurements.
Measure multiple distances	1. To measure three distances, tap 4 Point Dist .  2. To measure four distances, tap 5 Point Dist .  3. Position and anchor the beginning and ending markers.

General Doppler Measurement Methods

Use the configuration settings to customize the default measurement method for Doppler and the methods used for calculations.

The horizontal axis of the Doppler spectrum represents time and the vertical axis represents distance.

Prerequisite: The measurement function is active during Doppler.

To	Do This
Determine a pulsatility index from a system-generated trace	1. Tap PI Auto.  2. Position and anchor the first marker. 3. Position and anchor the second marker. The system draws a trace between the two markers, anchors a marker at the peak systole, and activates a marker for the end diastole. 4. Confirm the position at the end diastole and then anchor the marker.
Determine a pulsatility index from a user-generated trace	1. Tap PI Manual.  2. Position and anchor the first marker at the peak systole. 3. Position and anchor the second marker to trace one cycle on the spectrum. You draw a trace as you position the marker. When you anchor the marker, the system anchors a marker at the peak systole, and activates a marker for the end diastole. 4. Confirm the position at the end diastole and then anchor the marker.
Measure a velocity or frequency	1. Tap Velocity.  2. Position and anchor a marker.
Determine a resistive index	1. Tap RI.  2. Position and anchor the first marker on the peak systole. 3. Position and anchor the second marker on the end diastole.

To	Do This
Determine a flow volume using a diameter	This calculation requires a vessel diameter measurement in 2D-mode and a Doppler velocity trace. Perform the Doppler measurement from peak systole to the next peak systole or from the beginning of the cardiac cycle to the next end diastole. You can perform these measurements in any sequence. 1. Tap D-Flow.  2. Position and anchor the beginning and ending markers to measure the diameter of the vessel wall on the 2D-mode image. 3. Position and anchor the beginning and ending markers to include as many cardiac cycles as possible on the Doppler spectrum. Note: When you acquire the Doppler spectrum, extend the range of the Doppler gate to fit within the vessel of interest. Rotate ANGLE to adjust the Doppler angle parallel to the blood flow.
Determine a flow volume using an area	This calculation requires an area measurement in 2D-mode and a spectral Doppler trace. 1. Tap A-Flow.  2. Position and anchor the beginning marker of the ellipse to outline the inner lumen of the cross-sectional vessel on the 2D-mode image. 3. Rotate and size the first axis and then anchor the marker. 4. Adjust the length of the second axis to fit the ellipse to the structure and then anchor the marker. 5. Position and anchor the beginning and ending markers to include as many cardiac cycles as possible on the Doppler spectrum. Note: Perform the Doppler measurement from the peak systole to the next peak systole or from the beginning of the cardiac cycle to the next end diastole.
Determine the ratio of two velocity measurements	1. Tap Velocity Ratio.  2. Position and anchor the marker for the first velocity. 3. Position and anchor the marker for the second velocity.

To	Do This
Measure acceleration or deceleration	1. Tap Acceleration.  2. Position and anchor the beginning and ending markers along the horizontal and vertical axes. A dotted line indicates the slope.
Measure a time interval	1. Tap Time.  2. Position and anchor the beginning and ending markers.
Determine a velocity time integral trace	1. Tap VTI Trace.  2. Position and anchor the first marker. 3. Position and anchor the second marker to trace one cycle on the spectrum.
Determine a velocity time integral	1. Tap VTI.  2. Position and anchor the first marker. 3. Trace the structure. a. Position and anchor three or more additional markers. b. Press SET twice to anchor the ending marker.
Determine heart rate	1. Tap HR.  2. Position and anchor the first marker at the beginning of the cardiac cycle. 3. Position and anchor the second marker at the end of the cardiac cycle. Note: You can manually enter results for heart rate.
Determine a valve area by a pressure half-time	1. Tap PHT Guide.  2. Position and anchor the first marker on the spectrum at the peak velocity. 3. Position and anchor the second marker on the spectrum on deceleration of the mitral valve.

To	Do This
Determine a valve area on the irregular spectrum by pressure half-time	1. Tap PHT.  2. Position and anchor the first marker on the spectrum at the peak velocity. 3. Position and anchor the second marker on the spectrum on deceleration of the mitral valve.
Measure a pressure half time and deceleration time	1. Tap E/A Vel + DT. E/A 2. Position and anchor the first marker on the spectrum at the peak velocity. 3. Position and anchor the second marker on the spectrum on deceleration of the mitral valve.
Measure a rate of pressure change over time	1. Tap dP/dt. dP/dt 2. Position and anchor the beginning marker on the spectrum at the baseline of -100 cm/s. 3. Position and anchor the ending marker on the spectrum at the baseline of 300 cm/s.
Determine the ratio of two time measurements	1. Tap Time Ratio.  2. Position and anchor markers for two time measurements.

Follicle Measurements

(Available only for gynecology exams)

Use the configuration settings to select the measurement method for evaluating follicles and to enable or disable detection of follicles.

To perform a follicle measurement:

1. Activate the measurement function.
2. Tap the required measurement label.
3. Trace the structure using the method specified in the configuration settings.
4. Repeat steps 2 through 3 for each follicle.

To evaluate follicles with syngo Auto Follicle:

1. During a gynecology exam, acquire and freeze a 2D-mode image.
2. Tap **Auto Follicle** to activate the measurement function.

To	Do This
Perform a follicle measurement using Auto Follicle measurements	1. Tap the required measurement label on the touch screen. 2. Roll the trackball to position the measurement marker for each follicle and then press SET. The system traces the borders of the follicle and indicates one or two axes of the follicle depending on the Auto Follicle measurement method.  <i>Example of follicle measurement marker.</i> Note: The system displays a message if a follicle cannot be detected. You can reposition the measurement marker and repeat this step. 3. Press SET to confirm the measurement. The system updates the measured results for the measured follicle and activates a measurement marker for another follicle. 4. Repeat steps 2 and 3 for each follicle.
Edit a follicle measurement	Note: Before editing, press ESC if the marker is active. 1. Press UPDATE and then use the trackball to select the required follicle measurement for editing. 2. Position and then anchor the measurement marker for the selected follicle. The system displays the updated values in the measured results.
Deactivate follicle measurements	▪ Tap Auto Follicle again to turn off.

Performing Pediatric Hip Measurements

(Available only during 2D-mode imaging)

The pediatric hip measurement provides tools to detect dysplasia of the hip joint in infants by measuring the angles of the left and right hip and then determining the stressed and flexed percentages of femoral head coverage. The stability of the pediatric hip joint is determined by Graf classification using the alpha hip angle combined with the beta hip angle and acetabular coverage of the femoral head.

The hip angle measurement requires the placement of intersecting lines including a reference line to measure either the alpha angle only or both beta and alpha angles. Femoral head coverage is assessed by measuring the ratio of the acetabular depth to the diameter of the femoral head.

To determine a hip angle:

Prerequisite: Activate the measurement function and then select the measurement package for pediatric hip.

1. Tap a modifier for the measurement label.
 - **Right** indicates the right hip.
 - **Left** indicates the left hip.
 - **None** indicates a hip is not designated as left or right.
2. Tap the required measurement label.
 - **Rt Hip Angle** for the right hip angle
 - **Lt Hip Angle** for the left hip angle
 - **Hip Angle** for the hip not designated as left or right
3. Tap the required measurement tool.
 - **Hip (α)** determines only the alpha angle.
 - **Hip ($\alpha+\beta$)** determines the alpha angle and beta angle.
 - **Hip ($\beta+\alpha$)** determines the beta angle and alpha angle.

A reference line displays on the image screen. The mark on the line indicates the pivot point.

4. Position the reference line:
 - a. Position the line along the iliac bone.
 - b. Rotate **[Rotate]** to adjust the angle of the line.
 - c. Anchor the line.

5. Position the intersecting line to measure the alpha angle.
 - a. Position the line at the bony rim of the acetabulum.
 - b. Rotate **[Rotate]** to position the line along the lower limb of the iliac bone and intersect the reference line.
 - c. Anchor the line.
6. Position the intersecting line to measure the beta angle.

Note: The intersecting line to measure the beta angle is available only for **Hip ($\alpha+\beta$)** and **Hip ($\beta+\alpha$)**.

 - a. Position the line at the bony rim of the acetabulum.
 - b. Rotate **[Rotate]** to position the line along the labrum and intersect the reference line.
 - c. Anchor the line.

To determine a stressed or flexed percentage of femoral head coverage:

Prerequisite: Activate the measurement function and then select the measurement package for pediatric hip.

1. Tap a modifier for the measurement label.
 - **Right** indicates the right hip.
 - **Left** indicates the left hip.
 - **None** indicates a hip is not designated left or right.
2. Tap the required measurement label.
 - **Rt Flexed FHC** for the flexed percentage of femoral head coverage (right hip)
 - **Lt Flexed FHC** for the flexed percentage of femoral head coverage (left hip)
 - **Flexed FHC** for the flexed percentage of femoral head coverage (undesignated hip)
 - **Rt Stressed FHC** for the stressed percentage of femoral head coverage (right hip)
 - **Lt Stressed FHC** for the stressed percentage of femoral head coverage (left hip)
 - **Stressed FHC** for the stressed percentage of femoral head coverage (undesignated hip)
3. Tap the required measurement tool.
 - **%FHC** determines the percentage of femoral head coverage.
 - **d:D Ratio** determines the ratio of the acetabular depth to the diameter of the femoral head.

A reference line displays on the image screen. The mark on the line indicates the pivot point.

4. To determine the percentage of femoral head coverage using **%FHC**:
 - a. Position the reference line along the iliac bone.
 - b. Rotate **[Rotate]** to adjust the angle of the reference line and then anchor the line.
 - c. Position the circle to determine the femoral head diameter.
 - d. Rotate **[Size]** to adjust the size of the circle and then anchor the circle.
5. To determine the percentage of femoral head coverage using **d:D Ratio**:
 - a. Position and anchor the beginning and ending markers to measure the diameter of the femoral head.
 - b. Position and anchor the beginning and ending markers to measure the acetabular depth.

Calculated Results of the Pediatric Hip Tools

Results	Abbreviation
Alpha hip angle	α
Beta hip angle	β
Flexed percentage of femoral head coverage (undesignated hip)	Flexed FHC
Flexed percentage of femoral head coverage (right hip)	Rt Flexed FHC
Flexed percentage of femoral head coverage (left hip)	Lt Flexed FHC
Stressed percentage of femoral head coverage (undesignated hip)	Stressed FHC
Stressed percentage of femoral head coverage (right hip)	Rt Stressed FHC
Stressed percentage of femoral head coverage (left hip)	Lt Stressed FHC

Measuring Carotid Intima-Media Thickness with the Auto IMT Feature

(Not available during dual or live display format)

Auto IMT detects and measures the carotid intima-media thickness in a user-positioned region of interest on a frozen 2D image.

Use the configuration settings to show or hide the values in the measured results.

See also: For a list of measurement labels, refer to Appendix E in the System Reference.

To determine the carotid intima-media thickness:

Prerequisite: Select the carotid measurement package.

1. Freeze the image.
2. Press **CALIPER**.
3. Tap the required modifier and segment prefix modifier, for example, **Right** or **Prox**.
4. Tap the required measurement label.
The region of interest (ROI) displays on the 2D image. The dashed reference line is 1 cm from the ROI.
5. To preview temporary values for different locations of the vessel:
 - a. Position the ROI over the posterior wall of the common carotid artery and then press **UPDATE**.
The system traces the vessel and displays measured results for the selected ROI.
 - b. To adjust the IMT sensitivity, rotate **[IMT Sensitivity]**.
 - c. To exit the preview, roll the trackball.
6. Define the ROI.
 - a. Position the ROI over the posterior wall of the common carotid artery and then press a trackball key.
The system traces the vessel and displays two lines for measuring the thickness of the vessel.
 - b. To reposition the lines:
 - i. Press **SET** to edit the contour line.
The system indicates the editable line in green.
 - ii. Reposition the line and then press a trackball key to adjust the contour of the line.
7. Press **UPDATE** to complete the measurement.
The system confirms the contours, measures the anatomy, and displays the lines in yellow.

Calculated Results of Auto IMT Tools

Results	Abbreviation
Mean intima-media thickness	Mean
Maximum intima-media thickness	Max
Minimum intima-media thickness	Min
Intima-media thickness of the distal wall 10 mm proximal to the common carotid artery-carotid sinus junction	C10
The value is derived from the position of the region of interest as recommended by the Japanese Society of Ultrasound Medicine.	

Patient Reports

Patient reports arrange information into sections and include labeled measurement results and exam-specific information.

Use the configuration settings to customize the reports:

- Show or hide specific report sections and specify the sequence of the sections
- Provide site information, for example, a logo, boilerplate text, and a signature
- Provide login information
- Create predefined comments
- Customize lists of system-defined or user-defined anatomy structures and descriptions, and the sequence of the structures and descriptions in the lists

To access the modifiable report:

- Press **REPORT**.

Note: You can also view archived reports during review.

To	Do This
Display the exam specific information and measurement results	▪ Tap Results or click Exam Specific Information on the report menu. 
Display anatomy descriptions	▪ Click Assessment on the report menu. 
Display the biophysical profile	(Available only for obstetric exams) ▪ Click Biophysical Profile on the report menu. 
Display captures of my image screen captures and optional layouts	▪ Click Images on the report menu. 
Display the comments section	▪ Click Comments on the report menu. 
Open or close a report section	▪ Click the section name in the report.
Page through the report	▪ Click the scroll bar. Or, use the trackball to scroll upward or downward.
Display the image screen	▪ Press 2D.
Exit the report	▪ Press REPORT again. Or, tap Close.

Modifying Reports

To	Do This
Enter data in a text box	Click the text field and use the keyboard to enter the information.
Select an item from a list	Click the list and then click the item in the list.
Enter a date	- Click the calendar icon next to the date box. - Select the correct month and then click the date.
Edit a measurement	Note: If you edit a measurement required for a calculation, the system updates the calculated value. Click the measurement result and use the keyboard to replace the data. Edited measurements are marked with an asterisk (*). If an edited value is used in data averaging, the calculated average is marked with an asterisk.
Select the method for displaying a measurement value	- Select the required option from the drop-down list. Min uses the minimum value. Max uses the maximum value. Mean uses the mean of all measurement values. Last uses the last measurement value.
Insert an image	- Open the Images section. - Select the display layout from the drop-down list. - Click an image in the thumbnail panel. The system displays the image in the next available position.
Remove an image	Click X in the Images section.
Add a predefined comment	- Open the Comments section. - Click Load Comments. - Double-click the last row to add user-defined text to the list. - Enter a comment and then click Save.
Insert a predefined comment	- Open the Comments section. - Click Load Comments. - Double-click the required comment from the list. Or, click the required comment from the list and then click Add.
Edit a predefined comment	- Open the Comments section. - Click Load Comments. - Select the required comment and then click Edit. - Modify or delete the text. - Click Save. Or, press enter on the keyboard.

Tracing Areas of Interest

You can use the report drawing function to identify an area of interest on the body marker in the patient report.

Note: Use the prefix labels to indicate the right or left side of the anatomy.

To	Do This
Trace an area of interest	1. Open the Carotid Assessment or Prostate Assessment section in the report. 2. Scroll to the body marker in the report and then click Draw. The cursor becomes a pencil. 3. Use the trackball to trace the area of interest. 4. Press SET to stop tracing.
Remove a tracing	1. Click Delete. The cursor becomes an eraser. 2. Click a tracing to remove the tracing.
Remove all tracings	▪ Click Delete All.

Printing and Exporting Reports

To	Do This
View a report before printing	▪ Tap Print Preview. The system displays the preview of the report including registered patient information. Note: You cannot edit the report in print preview.
Redisplay the editable report	▪ Tap Print Preview.
Save the entire report as an image	▪ Tap Store Report.
Send the report to the configured printer	▪ Click Print Report.
Save the report as a portable document format (PDF)	1. Connect a USB storage device to a USB port. 2. Click Export PDF.

DICOM Structured Reporting


Note: DICOM structured reports are available only for the following exams: Cardiac, Obstetrics, Gynecology, and Vascular.

The ultrasound system generates a DICOM structured report (SR) for the current examination and stores the structured report to the study on the system's hard disk. You must send the structured report to a DICOM device to view the report.

Use the configuration settings to enable the system to send DICOM structured reports to DICOM devices.

Pediatric Hip Report

The pediatric hip report includes an assessment of the left and right hips and the Graf Sonometer, a depiction of the measured angles in a graph.

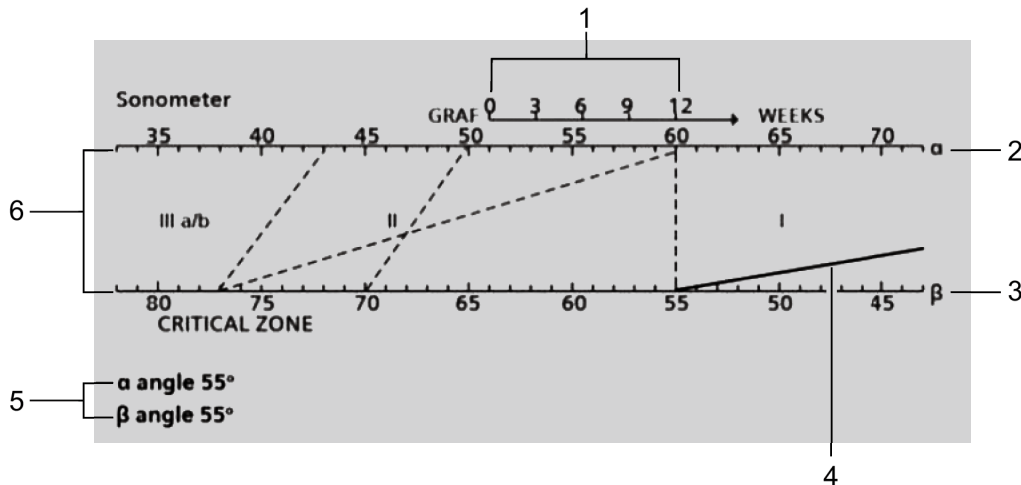

Example of the Graf Sonometer.

- 1 Acceptable normal range of alpha angle for patients younger than 12 weeks of age
- 2 Scale of degrees for the alpha angle
- 3 Scale of degrees for the beta angle
- 4 Solid line depicts the measured alpha and beta angles
- 5 Measured alpha and beta angles
- 6 Dashed lines indicate classification zones for the hip from dislocated to normal

To assess the left or right hip:

1. Press **REPORT**.
2. Click **Ped Hip Assessment**.
3. Click an option for assessment, for example, the osseous profile, osseous acetabulum, cartilaginous acetabulum, or an anatomical classification (type).
4. Select a classification from a drop-down list.

To view the measured angles in a graph:

(Available only when both the alpha and beta angles are measured)

1. Press **REPORT**.
2. Click the required report section:
 - **Right Hip Sonometer** for the right hip
 - **Left Hip Sonometer** for the left hip
 - **Hip Sonometer** for the hip not designated as left or right

C2 Obstetrical Measurements and Calculations

Overview.....	2
Obstetric Measurement and Calculation Labels.....	3
Controls for the Measurement Function	3
Making Obstetric Measurements.....	4
Measuring Multiple Fetuses.....	4
Performing Obstetric Measurements	5
Measuring Fetal Heart Rate.....	5
Using Auto OB.....	6
Determining an Amniotic Fluid Index	7
Obstetric Reports	8
Viewing Reports for Multiple Gestations.....	8
Modifying Measurements and Calculations in Reports.....	8
Modifying Descriptions in Reports	9
Viewing Doppler Results.....	9
Displaying Growth Analysis Graphs	10

Overview

The following obstetric exams are available on the ultrasound system:

- Exams for first trimester studies:
 - Early obstetric
 - Fetal echocardiography
- Exams for second and third trimester studies:
 - Obstetric
 - Advanced obstetric
 - Fetal echocardiography

Each exam type supports multiple gestations up to quadruplets.

You can generate growth analysis graphs to compare actual fetal growth to predicted growth patterns. You can create graphs to display measured results from a single exam, or you can include results acquired over a series of exams.

Use the configuration settings to customize the measurement function:

- Specify the system-defined measurement labels to display on the measurement menu, and customize the sequence, reference authors, and average or direct method for measurements and calculations
- Create user-defined measurement labels
- Create user-defined estimated fetal weight formulas and gestational age tables or formulas
- Create user-defined ratio formulas
- Create user-defined comments
- Customize the report

Obstetric Measurement and Calculation Labels

The obstetric exam types include measurement labels and calculation labels for 2D-mode, M-mode, and Doppler, and you can use the configuration settings to create user-defined measurement labels.


Obstetric calculations use tables or formulas designated by the reference authors specified in the configuration settings. If a measurement or calculation is outside the range specified by the selected reference, the system displays "xx w x d" instead of a result.

Calculation labels are not included on the measurement menu. The system performs the calculations when you make the required measurements. Calculation results are included with measurement results and in patient reports.

See also: For a list of measurement and calculation labels, refer to Appendix C in the System Reference.

Controls for the Measurement Function

You can access functions specific to the current operating mode or task using the selections in the left panel.

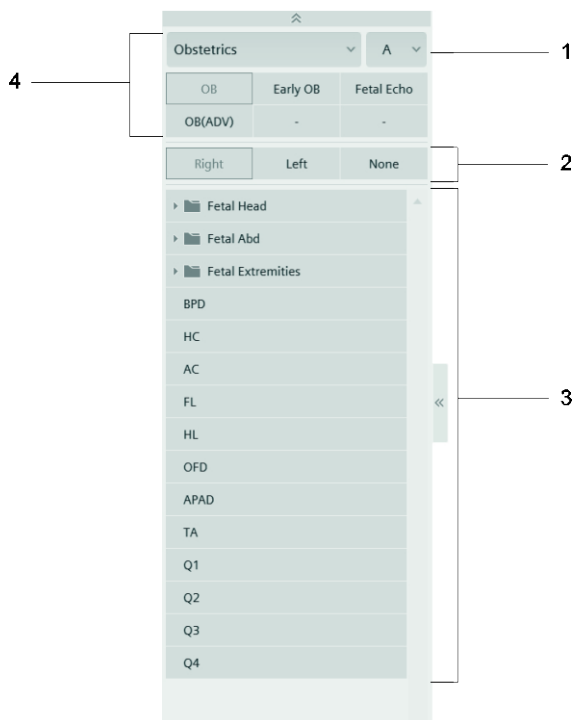

Example of measurement menu during an obstetric examination.

- 1 Identifier for multiple gestations
- 2 Prefix labels
- 3 Measurement labels

For obstetric parameters, the system displays the measurement result and the gestational age for the selected author.

- 4 Exam type

Displays an exam type. Click the arrow to select the required exam type from the list.

Making Obstetric Measurements

Prerequisite: The measurement function is active.

See also: For more information about general measurements, refer to Chapter C1 in this manual.

Measuring Multiple Fetuses

To	Do This
Measure multiple fetuses	1. Perform the measurements for the first fetus. The first measurements are assigned to fetus A. 2. Tap the multiple gestation indicator above the measurement menu and then select the next letter, for example, B, C, or D. 3. Perform the measurements for the next fetus. 4. Repeat steps 2 and 3 for each fetus.
Compare multiple fetuses	1. Perform the measurements for multiple fetuses. 2. Tap Fetus Compare. The system displays measured results for each fetus.
Display an estimated gestational age with crown rump length or biparietal diameter measurements	1. Measure a crown rump length or biparietal diameter. 2. Tap GA by CRL/BPD.
Assign measured results for APAD to TAD	1. Measure the anterior-posterior abdominal diameter (APAD) and the transabdominal trunk diameter (TAD). 2. Tap TAD/APAD.

Performing Obstetric Measurements

To	Do This
Measure the abdominal circumference	(Available only when AC with TAD/APAD Measurement is enabled in the configuration settings) 1. Tap AC. 2. Position and then anchor the first marker. The system displays an ellipse. 3. Position and then anchor the major axis. 4. Position and then anchor the second marker. 5. Roll the trackball to position the other axis. The system displays the measured results for the abdominal circumference (AC) with transverse abdominal diameter (TAD) and abdominal anteroposterior diameter (APAD). 6. Tap TAD/APAD to swap the values for TAD and APAD. 7. Press SET to confirm the measurement.
Measure the head circumference	(Available only when HC with BPD/OFD Measurement is enabled in the configuration settings) 1. Tap HC. 2. Position and then anchor the first marker. The system displays an ellipse. 3. Position and then anchor the major axis. 4. Position and then anchor the second marker. 5. Roll the trackball to position the other axis. The system displays the measured results for the head circumference (HC) with biparietal diameter (BPD) and occipital frontal diameter (OFD). 6. Tap BPD/OFD to swap the values for BPD and OFD. 7. Press SET to confirm the measurement.

Measuring Fetal Heart Rate

To measure fetal heart rate:

1. Tap **Fetal HR**.
2. Position and anchor the first marker at the beginning of the cardiac cycle.
3. Position and anchor the second marker at the end of the cardiac cycle.

Using Auto OB

The Auto OB measurement method positions measurement markers on detected fetal anatomy. You can reposition any marker, if necessary, and then anchor the final marker to determine the following measurements:

- Biparietal Diameter (BPD)
- Head Circumference (HC)
- Abdominal Circumference (AC)
- Femur Length (FL)
- Humerus Length (HL)

To	Do This
Activate Auto OB	1. Begin an OB exam in 2D-mode and then press FREEZE. 2. Press CALIPER to activate the measurement function. 3. Tap Auto OB. An orange dot indicates available Auto OB measurement labels. 4. Select the required Auto OB measurement label. The system displays the value in the measured results. 5. Reposition the measurement markers, if necessary. 6. Press SET to anchor the markers.
Deactivate Auto OB	▪ Tap Auto OB again.

Determining an Amniotic Fluid Index

The amniotic fluid index measurement determines a volume from four distance measurements, each measurement on a separate image in each of the four quadrants.

You can perform the procedure in 2D-mode or Dual mode. For 2D-mode, you must acquire a new image for each measurement. For Dual mode, you can make one measurement on each of the first two images and then acquire new images for the last two measurements.

To	Do This
Determine amniotic fluid index by measuring quadrants individually	1. Acquire and freeze a 2D image. 2. Tap a quadrant on the measurement menu: Q1, Q2, Q3, or Q4. You can make the quadrant measurements in any sequence. 3. Position and anchor the beginning and ending markers. The system displays the result for the selected quadrant and updates the calculated amniotic fluid index. 4. Repeat steps 2 and 3 until all quadrants are measured.
Determine the pocket depth of amniotic fluid	1. Acquire and freeze a 2D image. 2. Tap Amniotic Pocket on the measurement menu. 3. Position and anchor the beginning and ending markers to indicate the depth of the amniotic pocket.

Obstetric Reports

See also: For information on accessing, modifying, or printing a report, refer to Chapter C1 in this manual. Appendix C in the System Reference manual includes descriptions of sections in the report.

Viewing Reports for Multiple Gestations

The results for each fetus are included individually in the patient report.

To view a report for a specific fetus:

- Display the patient report and then select the indicator for the required fetus from the drop-down list.

Modifying Measurements and Calculations in Reports

When you modify either the clinical composite gestational age or the clinical estimated date of delivery, the system recalculates the date of last menstrual period.

Note: If you make modifications causing an estimated gestational age or estimated fetal weight to be outside the range specified by the selected reference, the system displays "xx w x d" instead of a result.

To modify a reference author in a report:

- Select a different author in the drop-down list for the required measurement.
The system updates the measurement and calculation results.
The system displays the last measured result only.

To select the method for displaying a measurement value:

- Select the required from the drop-down list.
 - **Min** uses the minimum value.
 - **Max** uses the maximum value.
 - **Mean** uses the mean of all measurement values.
 - **Last** uses the last measurement value.
 The system updates the measurement values and corresponding calculations in the report.

Modifying Descriptions in Reports

To assess fetal structures:

- Access the patient report and then click **Description** in left of the report page.

To assess fetal structures:

1. Access the patient report.
2. Tap **Assessment**.
3. Select the assessment from the drop-down list.

To	Do
Apply a predefined group of descriptions to all fetal structures in the report	▪ Click Load Defaults .
Clear the descriptions from all fetal structures in the report	▪ Click Clear All .
Define a new predefined group of descriptions	1. Select the required setting for each fetal structure. Note: You can include blank descriptions in the new predefined group. 2. Click Set as Default.

To create a biophysical profile:

1. Access the patient report.
2. Click **Biophysical Profile**.
3. Select a value of 0, 1, or 2 for each description:
 - Fetal Movement
 - Fetal Breathing Movement
 - Fetal Tone
 - Amniotic Fluid Volume

The system totals the four fields to calculate a score for the biophysical profile.

Viewing Doppler Results

(Available only for advanced obstetric exams)

You can view only Doppler measurements and calculations in the report or all measurements and calculations from all operating modes.

To include Doppler measurements and calculations in the report:

1. Access the patient report.
2. To display only measurements and calculations from Doppler, tap **Doppler Items**.
3. To display all measurements, including Doppler, tap **All Items**.

Displaying Growth Analysis Graphs

The graph section of the patient report includes growth analysis graphs. Growth analysis graphs display the upper limit, mean, and lower limit for predicted fetal growth patterns according to the selected reference for a measurement or calculation. The system plots the measurement results in the current exam on the vertical axis against the clinical age on the horizontal axis.

You can modify the labels and author references used to calculate the graphs. You can also combine results from previous exams with the current exam by manually entering measurement results or by linking to previous exams on the ultrasound system.

Use the configuration setting to create user-defined growth curves.

To view growth analysis graphs:

1. Tap **Growth Graph**.

Note: If the selection is not on the touch screen, access the report and scroll to view the growth graphs in the report.
2. To change the presentation of the graphs:
 - For advanced obstetric exams:
 - To view the graphs in a quad layout, tap **Display Quad**.
 - To view the graphs in a single layout, tap **Display Single**.
3. To include additional information in the graphs:
 - For early obstetric exams or obstetric exams:
 - To view the graph with data from previous studies, tap **Quik Trend**.
 - For advanced obstetric exams:
 - To view the graphs with the images, tap **Quick Grow Graph**.
 - To view the graphs with standard deviations, tap **Growing**.
 - To view the graphs with all standard deviations for estimated fetal weight, tap **Deviation All**.
 - To view the graphs with only a single deviation for estimated fetal weight, tap **Deviation Single**.
 - To view the graphs with an estimated gestational age, tap **USGA Data Graph**.
4. To modify the labels and references for the graph:
 - a. Select a different measurement, ratio, or calculation in the label drop-down list.
 - b. Select a different reference in the author drop-down list.

5. To link data from previous exams stored on the ultrasound system:
 - a. Click **Link Studies**.
 - b. Enter the patient name or patient ID and then click **Search**.
 - c. To exclude an exam, clear the check box next to the exam.
 - d. Click **Link**.

The system updates the graphs by combining the results from the previous exams with the current exam.
6. To manually enter results from previous exams:
 - a. Click **OB History**.
 - b. Click the calendar to enter the date of the exam.
 - c. Enter the measurement and calculation results.

The system does not check this data for accuracy or internal consistency.
 - d. Click **OK**.

C3 Cardiac Measurements and Calculations

Overview.....	2
Guided Measurements	2
Making Cardiac Measurements	3
Measuring a Heart Rate.....	3
Entering the Aliasing Velocity for Assessment of PISA	3
Assessing the Left Ventricular Function	4
Determining an End-Diastolic Volume or End-Systolic Volume	5
Using the eSie Measure Workflow Acceleration Package	6
Cardiac Reports.....	6

Overview

⚠ WARNING: Do not use the cardiac measurements and calculations package unless you are thoroughly trained in echocardiography and thoroughly familiar with the safe operation of the ultrasound system.

Imaging the heart necessitates a unique methodology for precision measurements and calculation of equations. 2D-mode, M-mode, and Doppler each have specific measurements and methods for cardiac exams. The system performs calculations when you complete the required measurements and assigns calculation labels to the results. Calculated results display in the Measured Results section of the screen and, along with labeled measurements, are transferred to the patient report.

Use the configuration settings to customize the measurement function:

- Select the default measurement method for each mode
- Select the default measurement method for each type of measurement
- For left-ventricle guided measurements, select measurement labels to include all the measurements for diastole and systole or to exclude specific labels from the guided measurement
- Customize the list of measurement labels for a non-guided measurement type
- Customize display options such as default measurement titles for each mode, the method used for data averaging, the method for determining Heart Rate, and the performing and referring physician information for the patient report
- Create user-defined cardiac measurement group and user-defined measurements label for each measurement group
- Derive the heart rate from a Doppler or M-mode signal or from the ECG trace
- Add eSie Measure measurements and series to user-defined workflow protocols

To select a cardiac measurement type:

1. Activate the measurement function.
2. Select the name of the measurement group from the measurement menu.
The corresponding measurement labels are in the measurement menu. The measurement methods for the active mode are available as soft key selections.

See also: For a list of cardiac measurement and calculation labels, refer to Appendix D in the System Reference.

Guided Measurements

Guided Measurements are system-defined sequences of individual measurements. You can perform measurements in the specified sequence or individually as non-guided measurements.

Making Cardiac Measurements

Prerequisite: The measurement function is active during system freeze.

See also: For more information about general measurements, refer to Chapter C1 in this manual.

Measuring a Heart Rate

You derive the heart rate from a Doppler or M-mode signal or from the ECG trace. In 2D-mode, you can manually edit the heart rate value.

To measure heart rate:

1. Tap **HR**.
2. Position and anchor the first marker at the beginning of the cardiac cycle.
3. Position and anchor the second marker at the end of the cardiac cycle.
4. To manually enter a heart rate, enter the heart rate in the dialog box.
The system displays the heart rate in yellow.

Entering the Aliasing Velocity for Assessment of PISA

You must manually enter the aliasing velocity value to properly assess mitral regurgitation, aortic regurgitation, and mitral stenosis of the proximal isovelocity surface area (PISA) during 2D-mode imaging and review as the ultrasound system does not display the aliasing velocity value acquired during live imaging.

To enter the aliasing velocity during 2D imaging or review:

1. Activate the measurement function.
2. Tap the **PISA (AR)**, **PISA (TR)**, or **PISA (MR)** measurement folder.
3. Tap the measurement label for the aliasing velocity.
4. Enter a value for the aliasing velocity and then click **OK**.
5. To exit without saving a value, click **Cancel**.

Assessing the Left Ventricular Function

Prerequisite: The measurement function is active during CINE playback of a 2D image or an M-mode sweep.

Analysis of the heart involves the assessment of the Left Ventricular (LV) function using 2D-mode or M-mode measurements.

You can select from the following volume formulas to calculate end-diastolic and end-systolic volumes.

- Formulas using distance and trace measurements:
 - Modified Simpson
 - Single Plane
 - Bi-Plane
 - Bullet
- Formulas using the disk volume method:
 - Simpson Single Plane
 - Simpson Bi-Plane
- Formulas using guided measurements for ventricular distances:
 - Cubed
 - Teichholz
 - Gibson

To perform a guided assessment:

1. Tap the required measurement group on the measurement menu.
2. Perform the measurements at end diastole:
 - a. Tap the diastole label on the measurement menu.
 - b. Position and anchor the marker on the posterior right ventricular wall.
 - c. Position and anchor the marker on the anterior interventricular septum.
 - d. During M-mode, position and anchor the marker on the posterior interventricular septum.
 - e. Position and anchor the marker on the endocardium of the left ventricular free wall.
 - f. Position and anchor the marker on the posterior left ventricular free wall.
3. Perform the measurements at end-systole:
 - a. Tap the systole label on the measurement menu.
 - b. Position and anchor the marker on the anterior interventricular septum.
 - c. Position and anchor the marker on the posterior interventricular septum.
 - d. Position and anchor the marker on the endocardium of the left ventricular free wall.
 - e. Position and anchor the marker on the posterior left ventricular free wall.
 - f. During M-mode, tap **LVET** and then position and anchor the measurement marker on the sweep.
4. Confirm the measured results.

To perform a non-guided assessment:

1. Tap **Cubed**, **Teichholz**, or **Gibson** on the measurement menu.
2. Tap any measurement label, except for the systole or diastole label.
3. Position and anchor the beginning and ending markers.
The system displays a checkmark next to the completed measurement and next to the diastole or systole series of measurements. The system highlights the next measurement in the group.
4. Continue to select measurement labels and perform measurements, as required.

Determining an End-Diastolic Volume or End-Systolic Volume

Prerequisite: The measurement function is active during CINE playback.

The Simpson Single Plane and Simpson Bi-Plane measurements use the disk volume measurement method.

To determine an end-diastolic volume or end-systolic volume:

1. Select a volume measurement label in the Measurement Menu.
2. Position and anchor the first marker to start a trace.
3. Trace the structure and anchor the ending marker.
If your ending marker does not meet the beginning marker, the system completes the trace and positions the axis line halfway between the two markers.
4. Position the axis line and anchor the ending marker.
5. Continue to select measurement labels and perform measurements, as required.

Using the eSie Measure Workflow Acceleration Package

(Not available for stored clips during review)

The eSie Measure workflow acceleration package reduces the keystrokes required to perform measurements. For example, the eSie Measure Teichholz activates the calipers and performs each required left ventricle measurement in diastole. When you confirm the measurements, the system copies the results to the patient report.

During M-mode and Doppler, the system displays measurements for at least three completed heart cycles or for all completed heart cycles if less than three measurements were performed.

To perform a measurement using the eSie Measure Function:

Prerequisite: An image or clip with at least one heart beat and an ECG trace.

1. Activate the measurement function.
2. Press **[eSie Measure]**.
3. Tap the required measurement group.
The system traces the image, displays a set of measurement results, and activates the last marker of the last measurement.
4. To manually perform an eSie Measure measurement or series, select the measurement label in the task pane and then position and anchor the markers.
5. If a message displays indicating that the system was unable to perform the measurement on the current image, click **OK**.
The system activates the measurement function.
6. To reposition a marker, select and then anchor the marker in the required location.

Cardiac Reports

See also: For information on accessing, modifying, or printing a report, refer to Chapter C1 in this manual. Appendix D in the System Reference manual includes descriptions of sections in the report.

D1 *syngo* Arterial Health Package

Overview.....	2
Terminology Used in This Chapter	2
Acquiring Clips for Use with the Software Program	3
Activating the Software Program	3
Reviewing Clips	4
Analyzing Clips	5
Modifying Settings	7
Accessing the CIMT Report	9
Distance Measurements in the Report	10
Clinical References.....	11

Overview

The *syngo* Arterial Health Package (AHP) software program provides a method of quantifying Carotid Intima-Media Thickness (CIMT). CIMT should be performed by a user who is experienced in CIMT measurement techniques.

Quantification of CIMT provides an estimate of vascular age for use in evaluating the potential risk of cardiovascular events in adults. The software program uses semi-automated border detection to determine the maximum and average thickness of the intima-medial layer of the carotid artery.

Terminology Used in This Chapter

Term	Definition
AHP	<i>syngo</i> Arterial Health Package (AHP) software program
CIMT	Carotid intima-media thickness
ASE	American Society of Echocardiography
ARIC	Atherosclerosis Risk in Communities study
CCA	Common carotid artery
CHD	Coronary heart disease
ECA	External carotid artery
ICA	Internal carotid artery
IMB	Intima-media boundary
IMT	Intima-media thickness
LBUL	Left bulb
LCCA	Left common carotid artery
LI	Lumen-intima
LIB	Lumen-intima boundary
LICA	Left internal carotid artery
MA	Media-adventitia
MAB	Media-adventitia boundary
ROI	Region of interest
RBUL	Right bulb
RCCA	Right common carotid artery
RICA	Right internal carotid artery

Acquiring Clips for Use with the Software Program

Note: Siemens Healthineers recommends following a description of the exam protocol, including acquisition of clips and images for analysis, from the American Society of Echocardiography (ASE) consensus statement.

Acquire clips from at least three angles of each blood vessel segment from both the left and the right sides of the body required to complete the selected protocol with an ECG active for adherence to the ASE consensus statement.

Activating the Software Program

You can activate the software program for a selected compatible clip in an active study or a previously saved study.

To activate *syngo* AHP:

1. Press **REVIEW**.
2. Select an image or clip.
3. Tap **AHP**.
The system displays an edit box for patient age and gender if there is missing patient information.
4. Review and analyze the clip.
5. To select another clip or image, double-click the required thumbnail.
6. To exit, tap **Exit**. Or, press **2D**.

Reviewing Clips

The thumbnail panel displays images and clips.


To review a clip:

- Double-click the required thumbnail.

To	Do This
Display the previous frame of the clip	Click Prior. Or, press the left arrow key on the keyboard. 
Start playback of the displayed clip	Click Run. 
Stop playback of the displayed clip	Click Stop. 
Display the next frame of the clip	Click Next. Or, press the right arrow key on the keyboard. 
Magnify the current clip or frame	Click Zoom +. Or, rotate ZOOM clockwise. 
Minimize the current clip or frame	Click Zoom -. Or, rotate ZOOM counterclockwise. 
Position the current clip or frame in the display area	Click Pan. 
Adjust the difference of the light and dark shades	Click Contrast. 
Lighten or darken shades	Click Brightness. 
Exclude a frame from quantification and analysis	- Click Toggle Exclude Frame.  Note: Excluding a frame does not delete the frame.  <i>Example of the symbol indicating the frame is unavailable for quantification and analysis.</i> - To include a frame, click Toggle Exclude Frame again.

Analyzing Clips

You can analyze any clip with only 2D-mode data and clips acquired with a linear transducer.

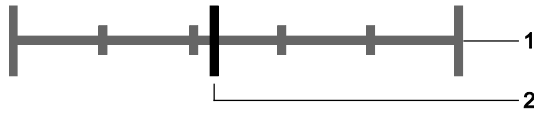

Example of displayed frame bar.

- 1 Frame bar (green line)
- 2 Reference position indicator (white line)

To analyze the contents of the clip:

1. Select a population group from the drop-down list on the top of the screen.
2. Select a vessel segment.
3. Select the orientation of vessel.

The software program outlines the CIMT region of interest (ROI) in red and displays a cursor on the left side of the ROI.

To	Do This
Position the CIMT ROI	▪ Click in the vessel near the far wall. Note: The upper edge of the ROI must be within the lumen of the vessel. The software program calculates and displays measurement values on the right side of the screen for the active frame.
Lengthen the CIMT ROI	1. Position the cursor on the outline of the CIMT ROI. 2. Press SET and then use the trackball to adjust the position of the CIMT region of interest. The system updates the measurement values.
Reposition the CIMT ROI without changing the length	▪ Press and hold CTRL and then press the right or left arrow key on the keyboard.
Edit the semi-automated trace borders	Note: Use the setup page to determine the ranges of color to apply to the "thickness" between the lumen-intima and media-adventitia or to only display lines to indicate the borders of the lumen-intima and media-adventitia. Prerequisite: Position the CIMT region of interest. 1. Click the required trace to edit and then click a key for the required edit: – Edit lumen-intima  – Edit media-adventitia  2. Position and anchor the markers along the boundary. A crossmark indicates the location of each edit point. 3. Position and anchor the last marker to conclude editing.  Example of the cursor during the edit function.

To	Do This
Discard the most recent edit	Click Undo Edit. 
Delete the outline of the CIMT ROI	Click Clear Borders. 
Trace a border manually	Note: Use the setup page to determine the ranges of color to apply to the "thickness" between the lumen-intima and media-adventitia or to only display lines to indicate the borders of the lumen-intima and media-adventitia. - Click a trace method: Manual trace lumen-intima  Manual trace media-adventitia  - Position and anchor the markers along the boundary. A crossmark indicates the location of each point. The software program connects the points to draw the line. Note: You must trace borders only drawn from left to right. - Click a measurement method. - Position and anchor the marker along the border of the vessel corresponding to the previous trace. - Position and anchor the ending marker after tracing the lumen-intima and media-adventitia.  <i>Example of the cursor during the manual trace function.</i>
Save CIMT data to the report	Click Save CIMT.
Return to the Carotid Intima-Media Thickness (CIMT) screen	Click CIMT Border. 
Magnify the selected region of interest (ROI)	Click Zoom to ROI. 

Modifying Settings

You can modify settings at any time. Use the setup page to assign the ranges of color to apply to the thickness between the lumen-intima and media-adventitia or to only display lines to indicate the borders of the lumen-intima and media-adventitia.

To modify settings for the region of interest:

1. Tap **Setup**.



2. Select the **Border** tab.

To	Do This
Determine the range of color to apply or to only display lines to indicate borders	1. Select the required option from the drop-down list in Display of CIMT area. CIMT applies color to the lines depicting the borders of the vessel boundaries or applies shades of color or a solid color to the "thickness" between the lumen-intima boundary and the media-adventitia boundary. A color bar displays the range of colors and a corresponding scale that indicates the average vessel thickness. You can assign the color to each color index if the shading option is selected. 2. Select a color scheme to apply to the selected CIMT area.
Specify the length or width of the CIMT region of interest	▪ Select a value in centimeters from the Search Length.
Specify the depth or height of the CIMT region of interest for the far wall	▪ Select a value in centimeters from the Far Wall Search Depth.
Visually indicate a segment thickness exceeds the user-defined threshold	▪ Select the Enable Maximum CIMT Threshold check box.
Specify the maximum threshold for the CIMT	▪ Select a value in millimeters from the Maximum CIMT Threshold. When reviewing a clip, the system displays a message indicating the CIMT threshold is exceeded and applies yellow highlights to the corresponding segments.
Enable the zoom feature for the region of interest	▪ Select the Zoom to ROI on Detection check box. Note: This feature is not applicable for detecting all frames.
Detect boundaries on all frames	▪ Select the Detect All Frames check box. Note: Detecting all frames may slow performance.
Display the CIMT region of interest	▪ Select the Show Search Region check box.
Omit the excluded frames from clip playback	▪ Select the Skip Excluded Frames on Playback check box.

3. To save changes, click **Save Page**.
4. To restore settings, click **Restore Page Defaults**.

To modify additional settings:1. Tap **Setup**.

To	Do This
Customize orientation descriptions for a protocol	1. Select the Define Protocol tab. The system displays a list of population sets and data values from the Atherosclerosis Risk in Communities (ARIC) study. 2. Ensure the corresponding protocol name for the specified vessel is displayed in the protocol box. 3. Click Copy Protocol. 4. Enter a name for the new protocol. The customized orientations are displayed on the measurements screen for the customized protocol.
Customize the protocol list	1. Select the Filter Protocols tab. 2. To display or exclude protocol types from the protocol list, select or clear the checkbox on the measurements screen.
Delete an existing protocol	(Available only with user-defined protocols) 1. Select the Define Protocol tab. 2. Select a protocol from the list. 3. Click Delete Protocol.
Set up a new report	1. Select the Report tab. 2. Click New in the Reporting Physician section. 3. Enter a name for a physician. 4. Click New in the Recommendations section. 5. Enter a title and description.
Remove a report setting	1. Select the Report tab. 2. Click Remove in the Reporting Physician or the Recommendations section.
Set up default units	1. Select the Units tab. 2. Select an option on the Units page and adjust the setting.

2. To save changes, click **Save Page**.3. To restore settings, click **Restore Page Defaults**.

Accessing the CIMT Report

Prerequisite: syngo AHP is active.

The report can display an overall vascular age and an overall risk for potential coronary heart disease. Risk analysis is based on information from the *Framingham Heart Study*.

See also: Clinical References, page D1-11.

The report can contain the following averages for each segment of the vessel:

- Maximum CIMT
- Mean CIMT
- Mean CIMT percentile (for age, sex and race)

The report can contain the following averages for both the right and left distal CCA segments:

- Mean CIMT
- Mean CIMT percentile
- Vascular age of the segment

The report can contain the following averages across all segments of the vessel:

- Composite mean CIMT
- Composite mean CIMT percentile
- Average vascular age

Note: To display the average vascular age and composite mean CIMT percentile values in the report, you must measure three angles of each segment. For ARIC protocols, you must measure one angle from each of the six segments.

- Framingham Model CHD Risk (%) (chronological age, blood pressure, cholesterol, diabetes, smoking)

To access the CIMT report:

- Tap **Reports**.



To	Do This
Edit or enter patient data relevant to <i>syngo</i> AHP calculations and Framingham assessment	1. Click Edit Patient Data. The system displays the screen with selections for entering risk factors. 2. Select the check box for Perform Framingham Assessment. 3. Enter the values for systolic and diastolic blood pressure. 4. Enter the values for both the total cholesterol and the HDL cholesterol levels. 5. Select the check box for additional risk factors.
Delete measurements	Note: You can delete measurements from the worksheet only. 1. Click Show Worksheet. 2. Click Edit Measurements. 3. Select the measurements to delete for a particular vessel, view, and wall. 4. Click Delete.
Enter sonographer initials, name of person preparing the report, reporting physician, recommendations, and a summary	1. Click Edit Summary. 2. Use the keyboard to enter your comments.
Export measurements or report data	▪ Click Export.
Display the report or worksheet	▪ To display the report, click Show Report. ▪ To display the worksheet, click Show Worksheet. ▪ To display the worksheet measurements screen, click Edit Measurements.

Distance Measurements in the Report

The software program provides the following measurements of the distances between the lumen-intima boundary and the media-adventitia boundary:

- Maximum CIMT distance for the currently displayed frame
- Mean CIMT distance for the currently displayed frame

Saved measurements display in the report and in the worksheet. You can delete the measurements only on the worksheet.

Clinical References

- Chambless LE, Heiss G, Folsom AR, Rosamond W, Szklo M, Sharrett AR, Clegg LX. "Association of Coronary Heart Disease Incidence with Carotid Arterial Wall Thickness and Major Risk Factors: The Atherosclerosis Risk in Communities (ARIC) Study, 1987-1993." *American Journal of Epidemiology* 146(6):483-494, 1997.
- Howard G, Sharrett A, Heiss G, Evans G, Chambless L, Riley W, Burke GL. "Carotid Artery Intimal-medial Thickness Distribution in General Populations as Evaluated by B-mode Ultrasound." *Stroke* 24(9):1297-1304, 1993.
- Stein JH. "Carotid Intima-Media Thickness and Vascular Age: You Are Only as Old as Your Arteries Look." *Journal of the American Society of Echocardiography* 17(6):686-689, 2004.
- Stein JH, Fraizer MC, Aeschlimann SE, Nelson-Worel J, McBride PE, Douglas PS. "Vascular Age: Integrating Carotid Intima-media Thickness Measurements with Global Coronary Risk Assessment." *Clinical Cardiology* 27(7):388-392, 2004.
- Stein JH, Korcarz CE, Hurst RT, Lonn E, Kendall CB, Mohler ER, Najjar SS, Rembold CM, Post WS. "Use of Carotid Ultrasound to Identify Subclinical Vascular Disease and Evaluate Cardiovascular Disease Risk: A Consensus Statement from the American Society of Echocardiography Carotid Intima-Media Thickness Task Force." *Journal of the American Society of Echocardiography* 21(2):93-111, 2008.
- Stein JH, Korcarz CE, Mays ME, Douglas PS, Palta M, Zhang H, LeCaire T, Paine D, Gustafson D, and Fan L. "A Semiautomated Ultrasound Border Detection Program That Facilitates Clinical Measurement of Ultrasound Carotid Intima-Media Thickness." *Journal of the American Society of Echocardiography* 18(3):244-251, 2005.
- Wilson PWF, D'Agostino RB, Levy D. "Prediction of Coronary Heart Disease Using Risk Factor Categories." *Circulation*, 97(18):1837-1847, 1998.


D2 eSie Left Heart

Overview of eSie Left Heart	2
Calculation Data and Measurements.....	3
Activating and Exiting eSie Left Heart	4
Working with eSie Left Heart	5

Overview of eSie Left Heart

⚠ WARNING: This application utilizes ultrasound data collected prior to analysis. The quality of the data collected and used as input could have an effect on the output of this application. Variability in ultrasound system performance, operator technique, patient characteristics, and other factors may affect the output. At all times, the clinician is advised to carefully review the output and confirm the information presented with clinical judgment and any other relevant sources of data.

eSie Left Heart is a clinical software program that detects the apical four-chamber or apical two-chamber view of the heart and then traces the endocardial borders (contours) of the left ventricle and left atrium. The program displays the contours on the end diastole (ED) and end systole (ES) images in a side-by-side format and generates measured results for end-diastolic volume, end-systolic volume, and ejection fraction. eSie Left Heart technology also includes methods for manually tracing the contours of the left ventricle and/or the left atrium.

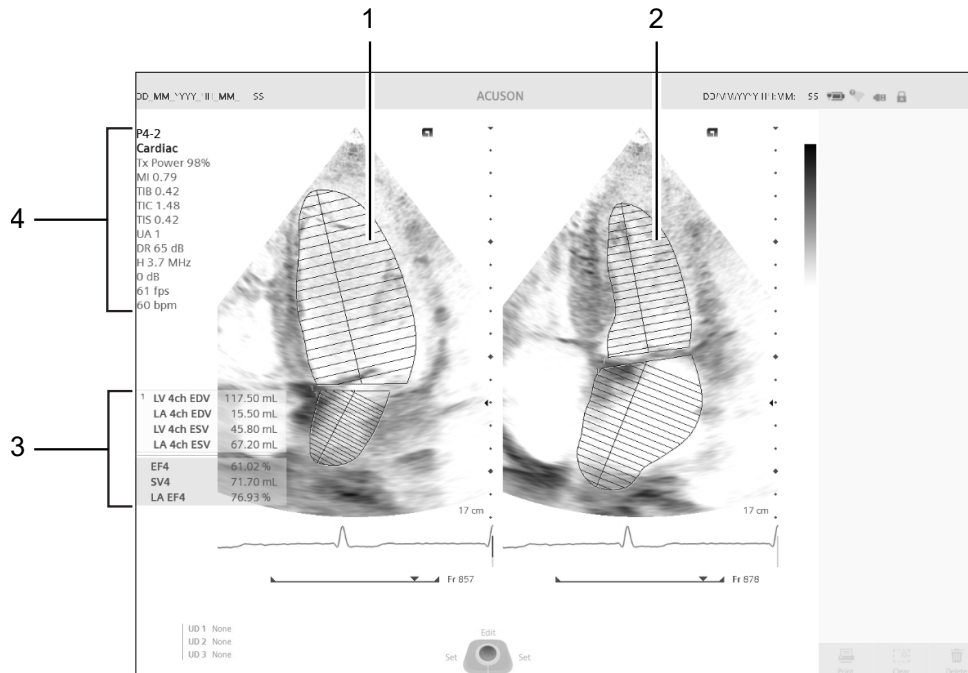

Example of eSie Left Heart screen.

- 1 Traced contours on the end diastole image
- 2 Traced contours on the end systole image
- 3 eSie Left Heart measured results
- 4 eSie Left Heart imaging parameters

Calculation Data and Measurements

When the traced contours for the ED and ES images are complete, eSie Left Heart generates the following measured results.

Note: Ejection Fraction Biplane calculation results are displayed when eSie Left Heart measurements have been completed for both the two-chamber and four-chamber views.

Measurement Label	Description
EF2	Left ventricle ejection fraction, apical 2-chamber view
EF4	Left ventricle ejection fraction, apical 4-chamber view
CO2	Left ventricle cardiac output for 2-chamber
CO4	Left ventricle cardiac output for 4-chamber
SV2	Left ventricle stroke volume for 2-chamber
SV4	Left ventricle stroke volume for 4-chamber
SI2	Stroke index for 2-chamber
SI4	Stroke index for 4-chamber
CI2	Cardiac index for 2-chamber
CI4	Cardiac index for 4-chamber
HR	Left ventricle heart rate
EDV (BP)	Left ventricle end-diastolic volume, bi-plane
ESV (BP)	Left ventricle end-systolic volume, bi-plane
LV 2ch EDV	Left ventricle end-diastolic volume, apical 2-chamber view
LV 2ch ESV	Left ventricle end-systolic volume, apical 2-chamber view
LV 4ch EDV	Left ventricle end-diastolic volume, apical 4-chamber view
LV 4ch ESV	Left ventricle end-systolic volume, apical 4-chamber view
EF (BP)	Left ventricle ejection fraction, biplane
LA EF2	Left atrium ejection fraction, apical 2-chamber view
LA EF4	Left atrium ejection fraction, apical 4-chamber view
LA SV2	Left atrium stroke volume for 2-chamber
LA SV4	Left atrium stroke volume for 4-chamber
LA EDV (BP)	Left atrium end-diastolic volume, bi-plane
LA ESV (BP)	Left atrium end-systolic volume, bi-plane
LA 2ch EDV	Left atrium end-diastolic volume, apical 2-chamber view
LA 2ch ESV	Left atrium end-systolic volume, apical 2-chamber view
LA 4ch EDV	Left atrium end-diastolic volume, apical 4-chamber view
LA 4ch ESV	Left atrium end-systolic volume, apical 4-chamber view
LA EF (BP)	Left atrium ejection fraction, biplane
LA Vol-I	Left atrium volume index

Activating and Exiting eSie Left Heart

You can activate eSie Left Heart on a real-time image containing only 2D-mode data, an ECG trace, and at least one heart cycle.

Note: eSie Left Heart supports clips captured in a reverse orientation using the left/right flip or up/down flip controls, or both.

To activate eSie Left Heart:

Prerequisite: An acquired and paused 2D clip in a full-screen display format.

1. Scroll to any frame in the required heart cycle, if necessary.
Note: If you select a frame in an incomplete heart cycle, eSie Left Heart uses the most recent complete heart cycle for analysis.
2. Tap **eSie LH**.
 The program displays the end diastole and end systole in a side-by-side display format and displays the measured results. The ED and ES images display the detected contours and major axes of the left ventricle and left atrium in yellow.
Note: The measured results display calculations only when the program detects the contours of the left atrium and the left ventricle or the contour of the left ventricle only. You can manually draw the contours or edit the program-generated contours.
3. To override the view label assigned by the system:
 - a. Tap **2CH** or **4CH**.
 - b. Tap the required measurement label for a two-chamber or four-chamber end-diastolic volume or end-systolic volume.
 The system reassigns the measured results and activates measurement labels in the task pane for the selected view.
4. To exit eSie Left Heart, tap **Exit**. Or, press **2D**.

Working with eSie Left Heart

eSie Left Heart displays contours only on the end diastole and end systole frames. You can manually adjust or delete and redraw the program-generated contours. Or, if the program could not detect the contours of the left ventricle or left atrium, you can manually draw the contours.

Note: Clicking the image, pressing **PRINT STORE**, or tapping **Edit** when the contours are displayed in green confirms the displayed the measured results. The editing button is activated on the touch screen only when you initially adjust and confirm the contours. When the contours are confirmed, the color turns to yellow.

To	Do This
Adjust the program-detected contours	1. Position the pointer on the required contour. The program displays the overall contour and segments of the contour in yellow. 2. To adjust a segment of the contour: a. For small adjustments, position the pointer closer to the contour and then click the image. b. For large adjustments, position the pointer farther from the contour and then click the image. The program displays the segment to be edited in green. 3. Reposition and anchor the segment using the trackball. 4. Repeat steps 2 and 3 to perform additional adjustments. When you move the pointer away from the contour, the program displays the contour in yellow.
Manually trace the contours of the left ventricle or left atrium	1. Tap the required location to manually trace the contours, for example 2Ch EDV, 2Ch ESV, 4Ch EDV, or 4Ch ESV. 2. Manually trace the contour by rolling the trackball. When you move the pointer away from the contour, the program displays the contours and major axes in green and displays the measured results on the screen.

To	Do This
Recalculate the measured results with a user-selected end diastole or end systole frame	1. Tap Previous frame or Next frame to view and select the required image (end diastole or end systole), if necessary. 2. Tap Set ED or Set ES. The program displays the detected contours on the end diastole and end systole frames and displays the measured results on the screen. 3. Tap ED Frame or ES Frame to view the user-confirmed end-diastolic or end-systolic frame.
Reposition the measured results on the screen	1. Click the measurement results. 2. Roll the trackball to the required location and then press SET.
Confirm the displayed contours and save the measured results to the patient report	▪ When the contours are displayed in yellow, confirm the contours by clicking the image, pressing PRINT STORE, or tapping Edit. The program displays the results of the method of disks on the displayed images.

D3 Transesophageal Transducer

Introduction to the Transesophageal Transducers	2
About the 5VT Transducer	3
Transducer Articulation Controls	4
Array Rotation Control	4
Flex Controls	5
Friction Brakes	7
Preparation for Use	9
Prior to First Use.....	10
Prior to Each Use	10
Inspecting the Transducer	11
Using a Transducer Sheath (Transducer Cover).....	12
Using Protective Drapes.....	12
Electrical Safety Considerations	12
Implantable Devices	12
Defibrillators	13
Punctures or Cracks	14
Electrosurgical Units.....	14
Exam Considerations	15
Bite Guards	15
Exam Procedures.....	15
Imaging with the 5VT Transducer	17
Changing the Scan Plane.....	18
Starting Position of the Transducer Array.....	19
Temperature Controls and Safeguards	20
Temperature Sensor	20
Managing the Lens Surface Temperature	21
Additional Thermal Safeguard	22
Caring for the Transesophageal Transducer	23
Cleaning and Storage.....	23
Approved Cleaners and Disinfectants for the Transesophageal Transducer	23
Cleaning and Disinfecting the Transducer.....	24
Cleaning a Transducer at the Point of Use.....	25
Cleaning a Transducer	26
Disinfecting a Transducer.....	27
Storing and Handling	29
Service and Repair	29
Preventing Injury and Equipment Damage.....	30
Transducer Technical Description	32

Introduction to the Transesophageal Transducers

⚠ WARNING: Do not use this transducer unless you are thoroughly trained in transesophageal echocardiography and are familiar with the orientation of cardiac images obtained through a TEE procedure. The transesophageal transducer should be used only by a licensed physician.

⚠ WARNING: Before attempting to use this transducer, you should be thoroughly familiar with the safe operation of the ultrasound system and this transducer. Refer to the content of this chapter and to the Safety and Care chapter in the Instructions for Use manual for ultrasound system and transducer-related safety information.


See also: For general ultrasound system and transducer-related safety information, refer to Chapter 2 and Chapter 6 in the Instructions for Use.

⚠ WARNING: The leakage current test for transesophageal transducers must be done prior to each use in order to reduce the likelihood of harm to the patient. Unsafe leakage current could result in electrical burns to the patient during use, including during defibrillation. Refer to the manufacturer's instructions included with the tester.

The transesophageal transducer can be introduced into the esophagus of the patient to obtain images of the heart structure. This procedure, known as transesophageal echocardiography (TEE), can provide virtually unobstructed views of the heart and surrounding tissue.

About the 5VT Transducer

The 5VT is a multi-frequency, multi-plane, phased sector array transducer.

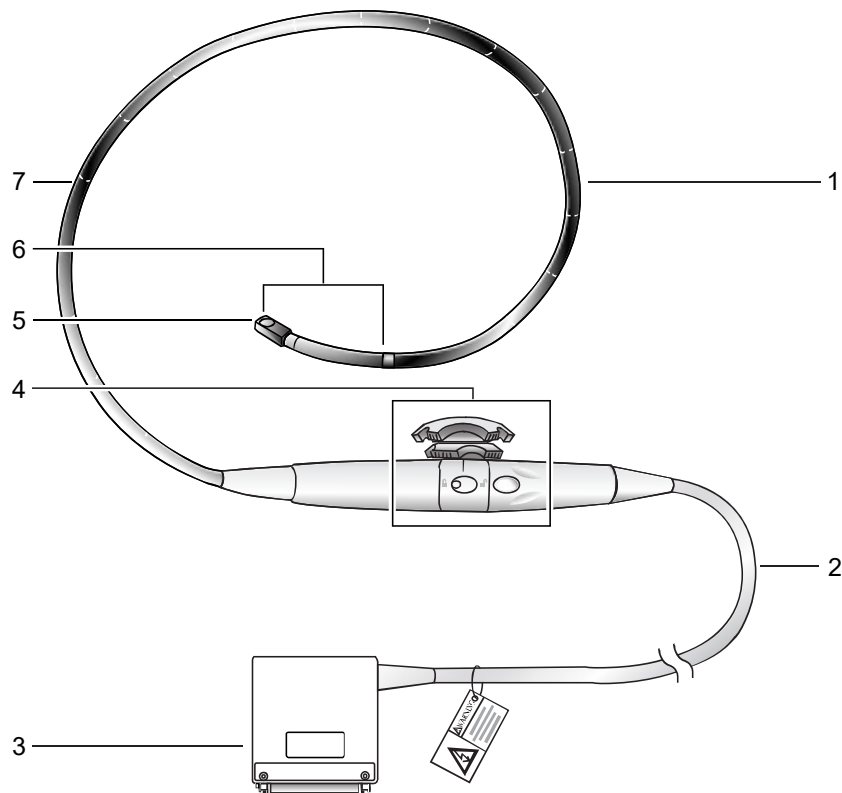


Controls located on the handle of the transducer allow anterior/posterior and left/right deflection of the distal tip and 180° rotation of the transducer array.

- 1 Flexible shaft
- 2 Cable
- 3 Transducer connector
- 4 Transducer controls
- 5 Distal tip (Contains rotating array)
- 6 Articulating section
- 7 Depth markings (in centimeters)

Transducer Articulation Controls

⚠ Caution: Do not use your hands or fingers to bend or manipulate the articulating section of the transesophageal transducer. To avoid severe damage to the articulating mechanism, use only the transducer controls to change the articulation angle.

The Array Rotation Control, Flex Controls, and Friction Brakes are located on the handle of the transducer.

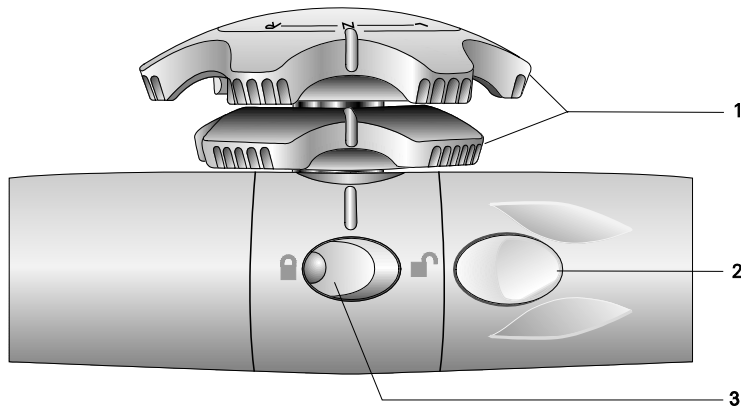

Example of transducer controls.


- 1 **Flex Controls**
Bend the articulating section of the transducer shaft
- 2 **Array Rotation Control**
Rotates the transducer array
- 3 **Friction Brakes** (one on each side of the transducer)
Locks the flex controls

Array Rotation Control

The array rotation control rotates the array within the distal tip of the transducer. The position of the array corresponds to the **scan plane**. The scan plane is adjustable from 0° to 180°.

Flex Controls

The flex controls manipulate the movement of the distal tip.

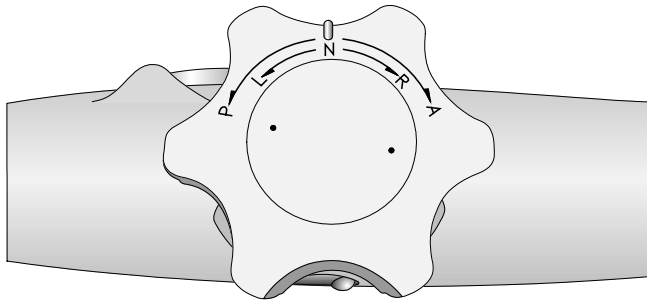

The Flex Controls are located on the handle of the transducer.

Anterior/Posterior Deflection

For anterior deflection (anteflexion) of the distal tip, rotate the anterior/posterior flex control toward the **A** marking to move the tip anteriorly.

For posterior deflection (retroflexion) of the distal tip, rotate the anterior/posterior flex control toward the **P** marking to move the tip posteriorly.


Right/Left Deflection

For right deflection of the distal tip, rotate the left/right flex control toward the **R** marking to move the tip right.

For left deflection of the distal tip, rotate the left/right flex control toward the **L** marking to move the tip left.

Neutral Position

To place the distal tip in the neutral position, rotate the anterior/posterior flex control and the right/left flex control to align the **N** marking and raised reference lines with the raised reference line on the handle of the transducer.

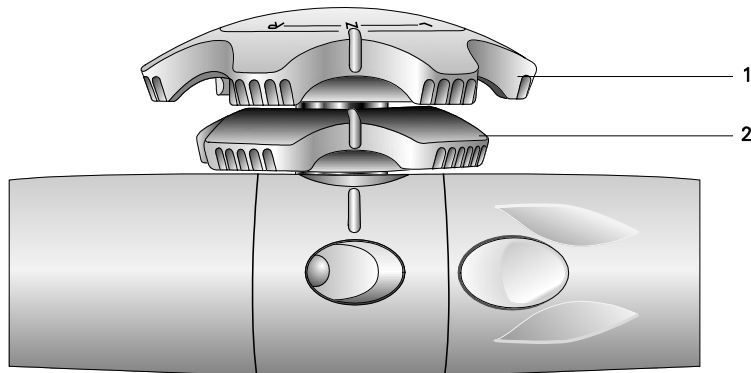


Flex controls in the neutral alignment.

- 1 Anterior/Posterior flex control
- 2 Left/Right flex control

Friction Brakes


The friction brakes lock the deflection positions of the distal tip. When the brakes are engaged, you can change the scan plane without losing acoustical contact.

A friction brake is located on each side of the transducer handle. The friction brakes are color-coded to match their corresponding flex control.

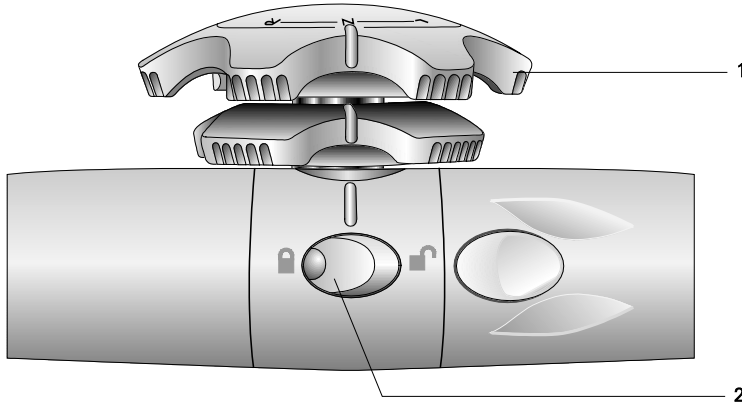

Example of transducer with locked/unlocked icons. Brake for Anterior/Posterior flex control is in the locked position.

- 1 Anterior/Posterior flex control (light gray)
- 2 Brake for Anterior/Posterior flex control (light gray)

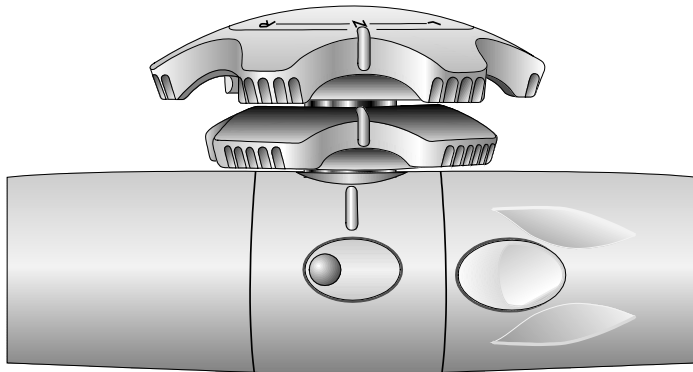


Example of transducer without locked/unlocked icons. Brake for Anterior/Posterior flex control is in the unlocked position.

To lock the anterior/posterior deflection position of the distal tip:

- Press the raised tactile indicator located on the button of the corresponding brake.


To unlock a brake:

- Press the flat portion on the button of the corresponding brake.

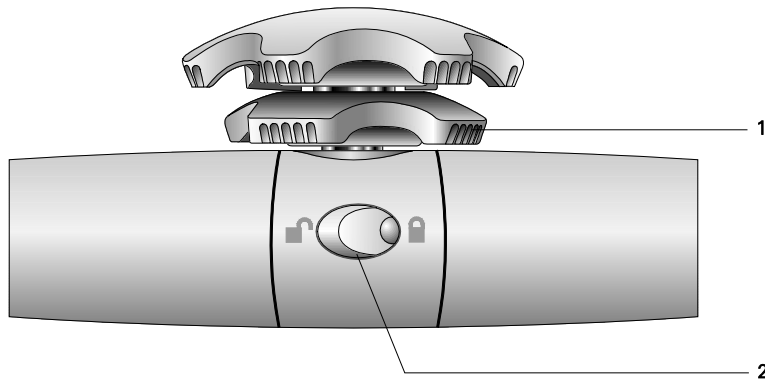

Example of transducer with locked/unlocked icons. Brake for Left/Right flex control is in the locked position.

- 1 Left/Right flex control (dark gray)
- 2 Brake for Left/Right flex control (dark gray)

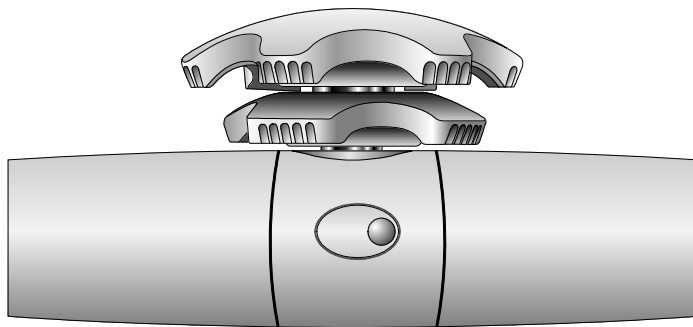

Example of transducer without locked/unlocked icons. Brake for Left/Right flex control is in the unlocked position.

To lock the left/right deflection position of the distal tip:

- Press the raised tactile indicator located on the button of the corresponding brake.

To unlock a brake:

- Press the flat portion on the button of the corresponding brake.

Preparation for Use

-  **WARNING:** When using an endocavity or intraoperative transducer with a BF or CF type applied part, the patient leakage currents may be additive.
-  **WARNING:** The leakage current test for transesophageal transducers must be done prior to each use in order to reduce the likelihood of harm to the patient. Refer to the manufacturer's instructions included with the tester.
-  **WARNING:** Prior to each use, inspect the endocavity or intraoperative transducer for signs of mechanical damage such as cracks, cuts, tears, perforations, or protrusions. Do not use the transducer if the transducer appears damaged in any way. A damaged transducer could cut the patient or user and compromise the electrical safety of the transducer, causing possible patient or user injury. Contact your local customer service representative.
-  **WARNING:** To prevent injury to you or the patient from electrical shock and burns when you connect or disconnect the transducer, do not touch the components on the transducer connector intended for insertion into the port on the ultrasound system. Grasp the transducer connector and avoid contact with the surface of the ultrasound system and any part of the transducer that touches the ultrasound system when disconnecting or connecting the transducer.

Prior to First Use

Proper precautions reduce risks to the patient and the equipment.

- Review the Safety and Care section.
- Thoroughly understand the features and operation of the transesophageal transducer, particularly the use of the flex controls and friction brakes.
- Clean and high-level disinfect the transesophageal transducer prior to its first use.
- Perform a leakage current test on the transducer. Refer to the manufacturer's instructions included with the tester.

Prior to Each Use

Review this checklist before each use of the transesophageal transducer to ensure patient safety, comfort, and confidence.

- Review patient history to ensure contraindications do not apply.
- Thoroughly inspect the transducer for mechanical damage.
- Familiarize yourself with the articulation controls.
- Ensure that the transducer has been properly cleaned and has been high-level disinfected, if disinfection is necessary.
- Perform a leakage current test on the transducer. Refer to the manufacturer's instructions included with the tester.
- Place a sheath on the transducer, if indicated.
- If necessary, place a drape over the transducer handle, cable, connector and the ultrasound system.

Inspecting the Transducer

To inspect the transducer before use:

1. Examine by look and touch the entire surface of the flexible shaft all the way to the distal tip for cuts, scrapes, protrusions, holes, dents, or cracks.
If you discover the surface of the transducer has been compromised, do not use the transducer. Contact your local customer service representative.
2. Test the array rotation control to ensure a smooth rotation.
3. Rotate the flex controls through the full range of motion and observe that the control movement is smooth and easy.
If the range of motion is tight or binding, or if the transducer controls make an unusual noise, do not use the transducer. Contact your local customer service representative.
If the distal tip appears to sag slightly when the flex controls are in the neutral position, this may be a sign of stretched or damaged control cables. Do not use the transducer. Contact your local customer service representative.
4. Test the brakes:
 - Engage the brakes (deflection locking mechanism) by pressing each ratchet button.
The brakes hold the distal tip at the appropriate articulation angle when locked.

Using a Transducer Sheath (Transducer Cover)

See also: For information about transducer sheaths, refer to Chapter 6 in the Instructions for Use.

⚠ WARNING: There have been reports of severe allergic reactions to medical devices containing latex (natural rubber). Health care professionals are advised to identify latex-sensitive patients and be prepared to treat allergic reactions promptly. For additional information in the U.S.A., refer to FDA Medical Alert MDA91-1.

For greatest patient and operator safety, cover a transducer with a transducer sheath. Use market-cleared transducer sheaths specifically designed for TEE applications. Follow the instructions provided by the manufacturer of the transducer sheath.

Disposal

While wearing protective gloves, remove the transducer sheath from the transducer. Dispose of the transducer sheath according to local, state, and regional laws and regulations for biohazardous waste.

Using Protective Drapes

Follow your hospital's guidelines regarding equipment use in the presence of infectious disease.

Electrical Safety Considerations

For safe use of this device, ensure an electrical safety procedure is in place for periodically inspecting the grounding system in the examination area and that this procedure is performed routinely.

Implantable Devices

⚠ WARNING: Ultrasound systems, like other medical equipment, use high-frequency electrical signals that can interfere with implantable devices such as pacemakers and implantable cardioverter-defibrillators (ICDs). When using a transesophageal transducer on a patient who has an implantable device, pay special attention to the performance of the implantable device. Stop the examination if there is interference.

Defibrillators

⚠ WARNING: Disconnect the transducer during defibrillation because a malfunction of the safety controls could otherwise result in electrical burns for the patient. Complete the leakage current test for transesophageal transducers in order to reduce the likelihood of harm to the patient during use.

The transesophageal transducer is not tested to meet certification requirements for a Type BF defibrillator-proof device per the IEC 60601-1 Standard for Safety of Medical Equipment.

Transesophageal transducers with this label are not certified as defibrillation proof:



Type BF Applied Part.

Examples of Risks

Follow the protocols of your clinic or hospital to determine the safest response for your patient before starting defibrillation: disconnecting the transesophageal transducer from the ultrasound system or removing the transesophageal transducer from the patient.

Transducer in the Patient during Defibrillation

- Defibrillation might not be as effective because the transesophageal transducer may absorb or divert the current required to restart the heart.
- If the electrical insulation surrounding the insertable portion of the transesophageal transducer becomes damaged during intubation, the patient could suffer burns to the esophagus. For example:
 - If the bite guard moves during intubation, the patient could bite into the inserted portion of the transesophageal transducer, exposing the transesophageal transducer's electrical circuitry to the patient's esophagus.
 - The inserted portion of the transesophageal transducer could become stressed during intubation, introducing cracks in the electrical insulation and exposing the transesophageal transducer's electrical circuitry to the patient's esophagus.

Transducer Removed from the Patient before Defibrillation

- Delay in defibrillation while the transesophageal transducer is removed from the patient.
- Injury to the esophagus, vocal cords, or teeth resulting from removing the transesophageal transducer too quickly.

Punctures or Cracks

Punctures or cracks in the outer layer of the transducer could expose the patient's esophagus to chassis leakage current and cause the patient to be at the chassis ground potential. If this occurs, the transducer will no longer meet the requirements of a Type BF classification. The transducer would no longer be a "floating" part and would fail a dielectric test. Provided that the internal grounding mechanism within the ultrasound system's power cord is intact and the cord is connected to a properly grounded wall outlet, there is little hazard to the user or patient from the system's leakage current. However, since the patient would now be at ground potential, leakage currents from other devices within the patient's environment could travel through the patient, compromising the safety to the patient or operator. If punctures or cracks in the outer layer of the transducer are observed, contact your local customer service representative.

Electrosurgical Units

⚠ WARNING: Use only isolated output electrosurgical units with the transesophageal transducer and disconnect the transducer when it is not in use. Failure to follow these procedures can result in esophageal burns for the patient, damage to the equipment, and unreliable data.

Use only electrosurgical units that have isolated outputs. If possible, use return-fault/ground-fault detection circuits, which provide extra protection.

Electrosurgical units and other operating room devices that do not have isolated outputs can introduce radio frequency electromagnetic fields or currents into the patient. The transducer is susceptible to these radio frequencies, which can interfere with the 2D-mode image and can completely override the color flow image, making it useless for diagnostic purposes.

In addition, any failure to an electrosurgical unit or other device, including the transesophageal transducer, could cause electrosurgical currents to return along the transducer's conductors. As a result, the patient could suffer from burns to the esophageal membranes. This arcing could also damage the transducer. To reduce the risk of leakage currents and electrosurgical interference, disconnect the transesophageal transducer from the ultrasound system when it is not in use.

Where isolation of electrosurgical units is in question, consult the manual for each electrosurgical unit or contact the biomedical engineering group.

Exam Considerations

Bite Guards

- ⚠ WARNING:** Follow the bite guard manufacturer's instructions to clean, disinfect, and rinse the bite guard if you use a reusable bite guard. Any disinfectant remaining on the reusable bite guard can burn the patient's mouth.

Bite guards are provided with the transducer for use with patients who will be awake during the exam. These guards are also available from third-party providers.

Do not use the provided bite guard with anesthetized patients. Use a separately purchased bite guard that has a space for taping an endotracheal tube to the guard.

Exam Procedures

- ⚠ WARNING:** Ensure that the friction brakes are unlocked and the flex controls are in the neutral alignment position when you insert or withdraw the transducer. Advancing or withdrawing the transducer while the transducer tip is bent and/or rigid can result in patient injury or death.
- ⚠ WARNING:** Do not use excessive force to advance or withdraw the transducer. Using excessive force to advance or withdraw the transducer can result in patient injury or death. If you encounter strong resistance during transducer articulation, discontinue the procedure. Identify and address the cause of the resistance before resuming the procedure. Withdraw and redirect the transducer as needed.
- ⚠ WARNING:** To reduce the risk of pressure necrosis, put the transducer in the neutral position when inserting or withdrawing it. Minimize pressure applied to the articulating section and distal tip. Do not let the distal tip displace tissue for more than five minutes.
- ⚠ WARNING:** Do not use the transesophageal transducer for any procedure requiring hyperextension of the neck. Improper patient positioning can cause paralysis of the vocal cord.
- ⚠ WARNING:** Always use a bite guard to prevent damage to the transducer from the patient's teeth, which could, as a result, create potential mechanical and electrical hazards for the patient. Using a bite guard also protects the patient's teeth.

Exam Procedures

Prerequisite: Connect the transducer to the ultrasound system and then activate the transducer prior to inserting the transducer into the patient.


1. Remove dentures, if present, and place them in a protected area until after the procedure.
2. Place a transducer sheath on the transducer, if needed.
3. Place a bite guard around the distal end of the transducer.
4. Ensure that both friction brakes are unlocked and both flex controls are in the neutral alignment position before inserting the transducer so that the transducer can bend to accommodate the path of the esophagus.
5. Place the transducer in the patient's mouth, and secure the bite guard into place before further advancing the transducer.
6. When the transducer has been inserted, press the array rotation control to change the scan plane.
7. As necessary, adjust the anterior/posterior and the left/right deflection using the flex controls. Position the transducer for optimal acoustical contact and the required viewing scan plane.
8. When the required scan plane is achieved and good acoustical contact is made, engage both the anterior/posterior and the left/right friction brakes.
9. Unlock both friction brakes and move both flex controls to the neutral position to acquire new echocardiographic views, which necessitate further advancement or withdrawal of the transducer, or change in the deflection.
10. Unlock both friction brakes and move both flex controls to the neutral position to remove the transducer.

Imaging with the 5VT Transducer

When you connect a transesophageal transducer to the system, the system displays a transducer-specific message.

Message	Displays When	Recommended Action
Refer to the warnings and precautions in the user manual before each use of the transesophageal transducer. Review instructions for leakage current testing and for inspection, proper use, and disinfection of the transesophageal transducer.	The transducer is selected.	Click OK .

During imaging with the 5VT transducer, information displays on the image screen regarding the rotation angle of the transducer array and regarding an indication of temperature at the lens surface of the transducer array.

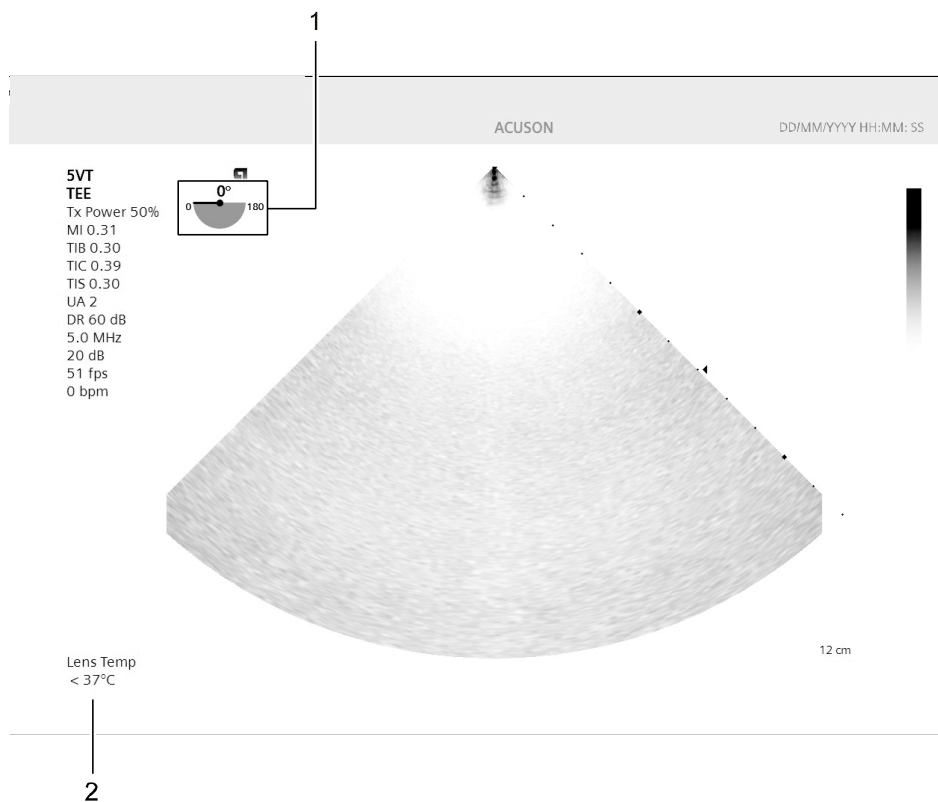

Example of the image screen during imaging with the 5VT transducer.


- 1 Scan plane icon indicating the approximate rotation angle of the scan plane
- 2 Temperature Indicator
Should normally display as **<37°C**.

Changing the Scan Plane

You can change the imaging plane of the transducer by using the array rotation control to rotate the transducer crystal within the distal tip of the transducer in a 0° to 180° range. The starting position for the transducer is 0° .

You can also change the imaging plane of the transducer using a selection assigned to the touch screen controls.

See also: For additional information about the touch screen controls, refer to Chapter 4 in the Instructions for Use.

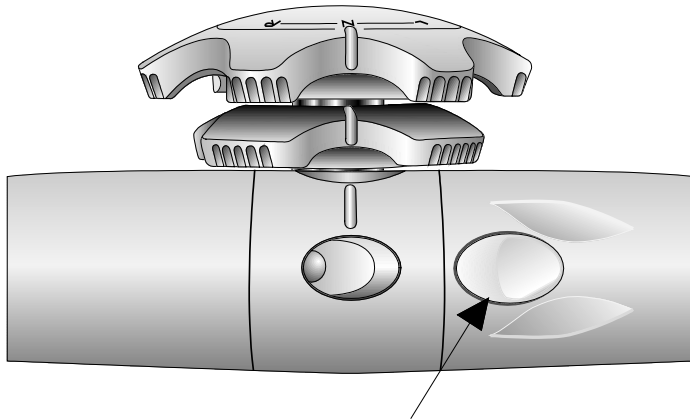

Example of the array rotation control is located on the handle of the transducer.


As the array orientation changes, the on-screen scan plane icon updates to indicate the relative direction of the array and the numerical value of the rotation angle.

The angle of 0° is equivalent to the transverse scan plane.	The angle of 90° is equivalent to the longitudinal scan plane.	The angle of 180° provides the mirrored transverse scan plane.

Examples of the scan plane icon. The array rotation angle is displayed numerically at the top of the icon.

To change the array orientation using the articulation control:

- Press the array rotation control to increase or decrease the array rotation angle. Press and hold the control to rapidly change the angle.

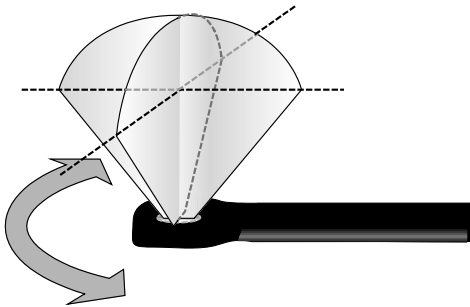

Example of array rotation on the transducer.

To change the array orientation in 90-degree increments:

- To rotate the image clockwise, tap the right arrow for **Rotate 90° (TEE)**.



- To rotate the image counterclockwise, tap the left arrow for **Rotate 90° (TEE)**.

**Starting Position of the Transducer Array**

When you first select the transducer, the system will automatically rotate the transducer array to the 0° position. If the system is unable to rotate the array to the 0° position, the system displays the following message:

The transducer rotation has not reset properly. Click OK to continue imaging at the current rotation. Or, to reset the transducer rotation, disconnect the transducer from the system, then reconnect it.

Choose a method:

- Click **OK** to continue imaging with the array in the current location.
The array rotation angle displays as --- at the top of the scan plane icon.
- Reselect or reconnect the transducer to reset the position of the array.

Temperature Controls and Safeguards

⚠ WARNING: Disconnect the transducer if the system displays the message indicating that the maximum thermal limit is exceeded. (The transducer temperature is 43°C or higher.) Failure to disconnect a transesophageal transducer that has exceeded the maximum thermal limit can burn the patient's esophagus. Reconnect the transducer after the transducer has cooled. Or, use a different transducer.

Because some concerns have been raised about possible thermal injury to the esophagus due to the local buildup of heat during transesophageal echocardiography, the transesophageal transducer incorporates a shut-off mechanism in the interest of patient safety in the unlikely event overheating should occur. For additional information on temperature limitation, refer to EN 60601-2-37 and IEC 60601-2-37.

Temperature Sensor

A temperature sensor mounted in the distal tip of the transducer monitors the temperature of the scanning array and displays the temperature on the screen.

When the temperature is 37°C or above, the system updates the value and color of the temperature indicator.

Temperature Indicator	Color of Temperature Indicator	Range of Lens Temperature
<37.0°C	Gray	Below 37.0°C
(Actual lens temperature)	Yellow	37.0°C to 40.9°C
(Actual lens temperature)	Orange	41.0°C to 42.9°C
43.0°C	Red	43.0°C
>43.0°C	Red	Above 43.0°C

Managing the Lens Surface Temperature

You can reduce the acoustic output to lower the lens surface temperature.

See also: Imaging Functions that Change Acoustic Output, Safety and Care, Chapter 2, Instructions for Use

To reduce the acoustic output:

- Lower the transmit power
- Change the scanning mode

Follow the instructions included in the on-screen messages if the lens surface temperature approaches, reaches, or exceeds the thermal limit.

Message	Displays when	Recommended Action
Transducer temperature is 42°C. Do you want to reduce the transmit power?	The lens surface temperature reaches 42°C , but the transmit power is not reduced.	▪ Choose a method: – Click Yes to automatically reduce the transmit power. – Click No to continue imaging while the lens surface temperature is 42°C. The system resets the override when you end the exam, select another transducer, or cycle power to the system. Note: You must disconnect the active transesophageal transducer before you can select a different transducer.
Transducer temperature is 42°C. Press OK to proceed with the current process. Or, wait for the transducer to cool.	The lens surface temperature reaches 42°C , but the transmit power is reduced.	▪ Choose a method: – Wait for the transducer to cool to 41°C. – Click OK to continue imaging while the lens surface temperature is 42°C.
Transducer temperature is 43°C. Imaging has been suspended. Allow transducer to cool to continue. Or, disconnect the transducer and reconnect after the transducer has cooled.	The lens surface temperature reaches 43°C (the maximum thermal limit).	When the lens surface temperature reaches 43°C, the system stops imaging and limits the available controls. 1. Disconnect the transducer. 2. Reconnect the transducer after the transducer has cooled. Or, use a different transducer. Note: The system is not operable until you disconnect the transesophageal transducer, or the transducer temperature cools down to 41°C.

Additional Thermal Safeguard

 **WARNING:** Disconnect the transducer if the system displays the message indicating that the acceptable thermal limit is exceeded. Remove the transducer from the patient. Failure to disconnect and remove a transesophageal transducer that has exceeded the maximum thermal limit can burn the patient's esophagus. Contact your customer service representative.

If the system's hardware detects the lens surface temperature is above the acceptable limit, the system stops imaging and displays the following message:

The transducer temperature is above the acceptable limit. Disconnect the transducer. Remove the transducer from the patient. Contact your Siemens Healthineers service representative.

You must observe the steps in the on-screen message.

Note: The system is not operable until you disconnect the transesophageal transducer or cycle power to the system.

Caring for the Transesophageal Transducer

See also: For general ultrasound system and transducer-related safety information, refer to Chapter 2 and Chapter 6 in the Instructions for Use.

Cleaning and Storage

You must properly clean and disinfect a transducer before the first use and after each use of a transducer.

Approved Cleaners and Disinfectants for the Transesophageal Transducer

- ⚠ WARNING:** The use of any cleaners and disinfectants other than those specified here may damage the transducer and, as a result, may create electrical hazards for patients and/or users.
- ⚠ WARNING:** Disinfectants and cleaning methods are listed for compatibility with product materials, not for biological effectiveness. Refer to disinfectant label instructions for guidance on disinfection efficacy and appropriate clinical uses.
- ⚠ WARNING:** Refer to the disinfectant manufacturer's instructions for guidance on proper rinsing of the disinfectant from the transesophageal transducer. Any disinfectant remaining on the transesophageal transducer can burn the patient's mouth or esophagus.

The Federal Republic of Germany prohibits the use of wipes for disinfection of transesophageal transducers and any transducer used in a semi-critical procedure.

See also: For an additional list of approved cleaners and disinfectants, refer to the Transducer Cleaners and Disinfectants Addendum. You can also access the most recent compatible cleaners and disinfectants on the Siemens Healthineers web site.

Transducer	Approved Cleaners	Approved Disinfectants
5VT	ENZOL	ANIOXYDE 1000
	Protex ULTRA Disinfectant Wipes	CIDEX OPA
	TEEZyme	gigasept FF (new)
	Transeptic Cleaning Solution	MetriCide OPA Plus ¹
		TD-5 ²
		TD-8 ²
		Tristel Trio Wipes System

¹For use as a high-level disinfectant, or for use with the ASTRA TEE reprocessing device according to the device manufacturer's instructions for use.

²For use with the TD-100 Automated TEE Probe Disinfector according to the device manufacturer's instructions for use.

Cleaning and Disinfecting the Transducer

-  **WARNING:** To avoid electrical shock and damage to the system, disconnect the transducer prior to cleaning or disinfecting.
-  **WARNING:** Never use Steam Heat or Ethylene Oxide for disinfecting the transesophageal transducer. These disinfection methods may damage the jacketing material, the transducer array, or the articulation controls, making the device unsafe to use on a patient.
-  **WARNING:** Never use iodine or solutions containing iodine for cleaning or disinfecting the transesophageal transducer. Iodine will cause degradation of the shaft material and of the distal tip resulting in a dangerous electrical hazard.
-  **Caution:** Never clean the transducer with bleach. Cleaning with bleach will damage the flexible shaft jacketing material, which will invalidate the manufacturer's warranty and could increase leakage current.
-  **Caution:** The distal tip is very sensitive to pressure. Care must be taken not to damage the tip while cleaning and disinfecting the distal tip.
-  **Caution:** Do not bend the flexible shaft into a curve of less than 30.5 cm in diameter.
-  **Caution:** The transducers have been designed and tested to be able to withstand high-level disinfection as recommended by the manufacturer of the disinfectant product. Carefully follow the disinfectant manufacturer's instructions.
-  **Caution:** Never immerse the entire transesophageal transducer in any solution. Immersion of the transducer handle could introduce fluid into the housing and cause damage to the controls or steering cables.
-  **Caution:** Do not clean the transducer with ethanol or methanol. Ethanol and methanol can damage the transducer.
-  **Caution:** Do not immerse the transducer above 100 cm from the distal tip; doing so may damage the controls.

Cleaning a Transducer at the Point of Use

Cleaning a transducer immediately after each use prevents soil and contaminants from drying on the surface of the transducer and facilitates subsequent cleaning.

To immediately clean a transesophageal transducer after each use:

1. Disconnect and remove the transducer from the ultrasound system.
2. Remove the transducer sheath.
3. Immediately clean the transducer in close proximity to the point of use.
 - a. Moisten a clean gauze pad or lint-free cloth with water or an approved cleaner. Or, use an approved cleaning wipe or sponge.
See also: For a list of approved cleaners, refer to page D3-23.
 - b. Carefully wipe the transducer surface intended for contact with the patient. Avoid touching the electrical components on the connector.
4. Perform subsequent cleaning of the transducer according to the next procedure.

Cleaning a Transducer

Cleaning a transducer removes soil and contaminants from the surface of the transducer and prepares the transducer for disinfection.


To clean a transesophageal transducer:

Prerequisite: The transducer is disconnected and removed from the ultrasound system. The transducer sheath is removed from the transducer.

1. Select an approved cleaner. Follow the manufacturer's instructions for concentration and temperature of the solution and contact time when cleaning the surface of the transducer.
See also: For a list of approved cleaners, refer to page D3-23.
2. Clean the distal tip and flexible shaft of the transducer using one of the following cleaning methods. Avoid touching the electrical components on the connector.
 - If immersion is required, submerge the distal tip and flexible shaft up to the 100-cm mark. Do not allow the handle, controls, cable, or connector to become wet.
 - If wiping is required, carefully wipe the distal tip and flexible shaft.
 - If spraying is required, do not spray the distal tip and flexible shaft near the ultrasound system. Never spray the transducer connector.
3. Clean the handle, controls, cable, and connector of the transducer. Keep the electrical components on the connector dry.
 - a. Moisten a clean gauze pad or lint-free cloth with an approved cleaner. Or, use an approved cleaning wipe.
 - b. Carefully wipe the handle, controls, cable, and connector.
4. Rinse the distal tip and flexible shaft up to the 100-cm mark according to the cleaner manufacturer's instructions. Do not rinse the handle, controls, cable, and connector. Keep the electrical components on the connector dry.
Note: If the cleaner manufacturer does not provide rinsing instructions, thoroughly rinse the distal tip and flexible shaft up to the 100-cm mark with deionized water at room temperature.
5. Disinfect the transducer.

Disinfecting a Transducer

Immediately after cleaning, disinfect a transducer to destroy pathogens and other microorganisms and to prepare the transducer for storage or the next clinical use.

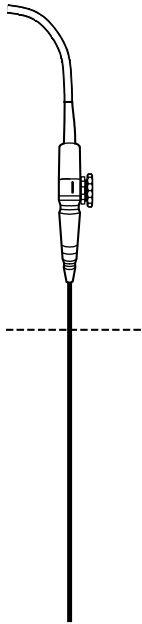

Submerge only to the 100-cm mark.

To disinfect a transesophageal transducer:

Prerequisite: Clean the transducer.

1. Select an approved disinfectant for high-level disinfection. Follow the disinfectant manufacturer's instructions for concentration and temperature of the solution and contact time when high-level disinfecting the surface of the transducer.
See also: For a list of approved disinfectants, refer to page D3-23.
2. Disinfect the distal tip and flexible shaft of the transducer using one of the following disinfection methods. Avoid touching the electrical components on the connector.
 - If immersion is required, submerge the distal tip and flexible shaft up to the 100-cm mark. Do not allow the handle, controls, cable, or connector to become wet.
 - If wiping is required, carefully wipe the distal tip and flexible shaft.
 - If spraying is required, do not spray the distal tip and flexible shaft near the ultrasound system. Never spray the transducer connector.
 - If a reprocessing device is required, use a device compliant with ISO 15883 and follow the reprocessing device manufacturer's instructions.

3. Disinfect the handle, controls, cable, and connector of the transducer. Keep the electrical components on the connector dry.
 - a. Moisten a clean gauze pad or lint-free cloth with an approved high-level disinfectant. Or, use an approved disinfectant wipe.
 - b. Carefully wipe the handle, controls, cable, and connector.
4. Rinse the distal tip and flexible shaft of the transducer up to the 100-cm mark according to the disinfectant manufacturer's instructions. Do not rinse the handle, controls, cable, and connector. Keep the electrical components on the connector dry.

Note: If the disinfectant manufacturer does not provide rinsing instructions, thoroughly rinse the distal tip and flexible shaft up to the 100-cm mark with deionized water at room temperature.
5. Follow the disinfectant manufacturer's instructions for drying the transducer. Do not use heated air.

Note: If the disinfectant manufacturer does not provide drying instructions, dry the transducer with a clean, lint-free, soft cloth.
6. Apply an approved disposable tip protector to the distal tip before storage.
7. Store the transducer in a clean, dry, protected environment. The transesophageal transducer is fragile and should not be stored in the transducer holders on the system.

Storing and Handling

The transesophageal transducer is a mechanically sensitive device. The distal tip is especially sensitive to bumps, scrapes, and jarring. Take extreme care when handling and storing the transducer. Never move the distal tip by hand. Always use the flex controls.

When you store the transesophageal transducer, ensure that the friction brakes are off, the distal tip is straight, and the cables and flexible shaft are free of kinks.

Use an approved disposable tip protector on the distal tip during storage.

Store the transducer in a clean, dry, protected environment. The transesophageal transducer is fragile and should not be stored in the transducer holders on the system.

Transporting the Transducer

Use the case for shipping or for moving the disinfected transducer from one site to another.

Take care that the distal tip is straight and that the cable is free of kinks as you carefully place the transducer in the foam padding. Before you close the lid, ensure that nothing protrudes from the case.

Service and Repair

Do not attempt to repair or alter any part of the transesophageal transducer. Contact your local customer service representative immediately if your transducer appears to be damaged or malfunctions in any way.

Preventing Injury and Equipment Damage

Following the preventive measures outlined below will help ensure patient safety.

Preventive Measure	Problem	Potential Harm to Patient/Potential Damage
Place the transesophageal transducer's flex controls in the neutral alignment position, without applying the friction brakes, when you insert or withdraw the transducer. Do not force the transducer.	Improper insertion or withdrawal	Esophageal cuts, bruising, bleeding, ligament damage, perforations
Inspect the transesophageal transducer prior to each use to ensure that it is functioning properly and that it has no cuts, tears, fraying, loose parts, or rough edges.	Mechanical damage, electrical damage, leakage current	Severe trauma, esophageal cuts, bleeding, perforation, electrical burns, serious electrical hazards
Use only isolated output electrosurgical units with the transesophageal transducer. Disconnect the transducer when it is not in use.	Non-isolated electrosurgical units	Electrical burns
Adjust the Transmit Voltage (dB or %) as low as clinically useful. Observe the on-screen temperature warning messages and respond accordingly.	Improper transducer temperature and/or acoustic output	Esophageal burns
For maximum safety, disconnect the transesophageal transducer from the system prior to defibrillation rather than relying on the transducer's built-in safeguards.	Defibrillation	Electrical burns
Maintain a regimen of electrical safety checks.	Electrical damage, leakage current	Esophageal burns, electrical damage
Put the transducer in the neutral position, without either friction brake applied, when inserting or withdrawing it. Minimize pressure applied to the articulating section and distal tip. Do not let the distal tip displace tissue for more than five minutes.	Pressure necrosis	Permanent damage to esophageal lining

Preventive Measure	Problem	Potential Harm to Patient/Potential Damage
Do not use the transesophageal transducer if a patient has esophageal varices or masses or strictures are present.	Esophageal Varices	Excessive bleeding
Do not use the transesophageal transducer for any procedure requiring hyperextension of the neck.	Improper patient position	Paralysis of the vocal cord
Be aware that a patient with a loose tooth could lose the tooth in a TEE procedure.	Improper assessment of dental risk	Tooth loss
Use a bite guard.	No dental protection, leakage current	Dental damage, electrical damage
Follow the bite guard manufacturer's instructions to clean, disinfect and rinse the bite guard if you use a reusable bite guard. Refer to the disinfectant manufacturer's instructions for guidance on proper rinsing of the disinfectant from the transesophageal transducer	Disinfectant residue on the bite guard and transesophageal transducer	Burns in the mouth and esophagus
Never use Steam Heat or Ethylene Oxide for disinfecting the transesophageal transducer. Do not use any disinfectants other than those specified. Do not use iodine or solutions containing iodine for cleaning or disinfecting the transesophageal transducer.	Improper cleaning resulting in electrical damage	Esophageal burns
Immerse the transesophageal transducer no longer than the time required by the manufacturer for high-level disinfection. Never immerse transducer controls or the system connector in liquids.	Improper disinfection, exposing electrical workings to liquids	Nosocomial infection, electrical damage affecting image quality, electrical safety, or mechanical operation
Do not force distal tip articulation. Always use the controls to alter articulation.	Forcing the deflection control	Stretching or breaking deflection controls
When the transesophageal transducer is not in use, store the transducer in a clean, dry, protected environment. Do not store the transducer in a tight coil.	Improper storage	Damage to array elements, cuts and abrasions to the flexible shaft

Transducer Technical Description

Transducer Type:	Phased Array, multiplane
Elements:	64
Array Rotation Angle:	0° – 180°
Distal Tip Articulation:	120° anterior, 90° posterior 45° right, 45° left
Maximum insertion length:	1080 mm
Maximum insertion width:	15 mm
2D-mode Frequencies:	3.6 MHz, 4.0 MHz, 5.0 MHz, 5.7 MHz
M-mode Frequencies:	3.6 MHz, 4.0 MHz, 5.0 MHz, 5.7 MHz
Color Doppler Frequencies:	3.1 MHz, 3.8 MHz
PW Doppler Frequencies:	3.1 MHz, 3.8 MHz
CW Doppler Frequencies:	2.9 MHz, 3.3 MHz
Imaging Modes:	2D-mode Color Doppler (color, power, DTI, DTE) M-mode Color M-mode (color) Pulsed Wave (PW) Doppler Pulsed Wave (PW) Doppler Tissue Imaging (DTI) Steerable Continuous Wave (SCW) Doppler
Maximum Image Field of View:	90°
Image:	Sector
Operating Environment Temperature: Relative Humidity:	 10°C to +40°C 10% to 80%
Storage Environment Temperature: Relative Humidity:	 -10°C to +50°C 10% to 85%
Degree of protection against electrical shock:	 Transducers with this label are certified as Type BF (not defibrillator proof) per IEC 60601-1 Standard for Safety of Medical Equipment
Degree of protection against ingress of fluid:	IPX8
Measurement Accuracy Rotation angle of the transducer array: Temperature: Accuracy:	 $\leq \pm 5^\circ$ 35°C to 47°C $< \pm 1^\circ\text{C}$

D4 *syngo* Velocity Vector Imaging


Overview.....	2
Selecting a Clip for Analysis	3
Analyzing a View	4
Analyzing a Different Clip in the Same View	6
Analyzing a Clip in a Different View.....	6
Setting Up a Clip for Processing	7
Displaying M-mode Data	8
Drawing and Modifying a Contour	9
Modifying a Contour	12
Processing a Contour and Viewing Results.....	13
Viewing Results in a Bullseye Plot.....	15
Analysis Page: Contour	16
Editing and Deleting a Trace	16
Analysis Page: Settings	18
Analysis Page: Display Buttons	20
Analysis Page.....	21
M-Mode	21
Contour.....	21
Segmental Analysis.....	21
Segmental Analysis Results	22
Viewing Segmental Analysis Results.....	23
Analysis Display Information.....	24
Saving Analysis Results to a Report	26
Viewing Analysis Results.....	26
Results by Parameter.....	26
Exporting Images, Clips, and Data.....	27

Overview

⚠ WARNING: *syngo* Velocity Vector Imaging technology (*syngo* VVI) utilizes ultrasound data collected prior to analysis. The quality of the data collected and used as input could have an effect on the output of this program. Variability in ultrasound system performance, operator technique, patient characteristics, and other factors may affect the output. At all times, the clinician is advised to carefully review the output and confirm the information presented with clinical judgment and any other relevant sources of data.

syngo VVI is a clinical software program used to visualize, measure, and assess myocardial motion and mechanics using acquired 2D clips. The program tracks and estimates tissue velocity and other motion and deformation parameters at selected points on a user-defined outline (contour) of a structure.

This program illustrates tissue velocity by generating velocity vectors (arrows) on the contour. These velocity vectors are the source of the derived quantitative data, such as velocity and strain rate. You can use this program to analyze rotation, displacement, and radial strain of the left ventricle, right ventricle, and left atrium in pediatric and adult studies.

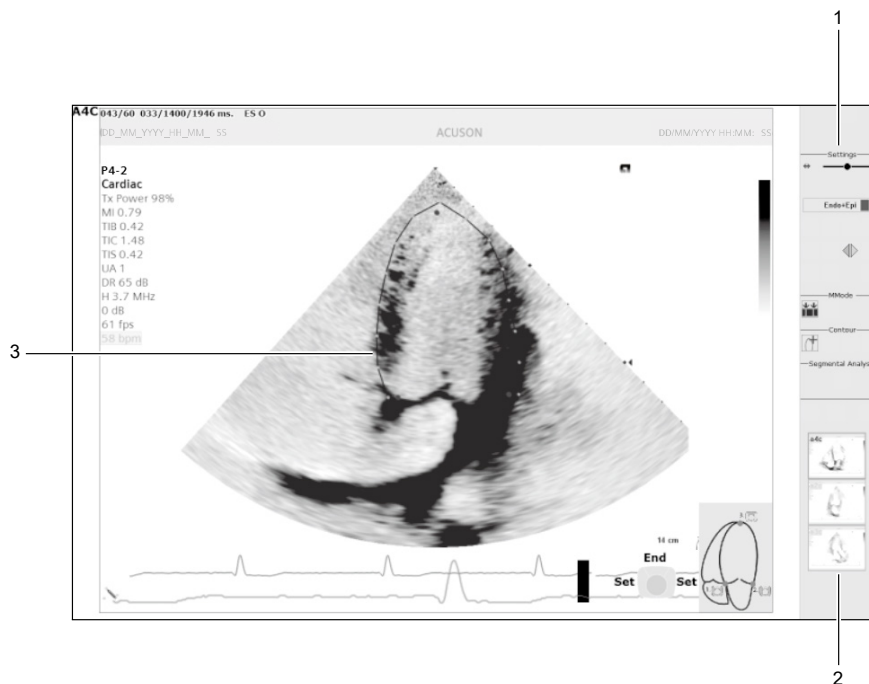

Example of an anatomy boundary definition on the Contour screen.

- 1 Selections for displaying and editing the contour
- 2 Thumbnail of a clip selected for analysis
- 3 Contour

Selecting a Clip for Analysis

Note: You can activate *syngo* VVI only if you select a valid clip. The system can load up to 10 clips.

You can activate this program using a 2D-mode clip that contains at least six frames of 2D-mode data. If the clip contains more than 700 frames, the program loads the first 700 frames acquired.


This program also supports the following functionality:

- Playing back and editing clips using the control panel
- Capturing static images

You cannot load 2D color clips, secondary capture images or clips, M-mode or Doppler images, static captures, 2D color dual, 2D live dual clips, 2D dual updated dual clips, and 2D clips with the trigger function enabled.

If you load clips with annotations, you can view the annotations, but you cannot activate the annotation feature.

Analyzing a View

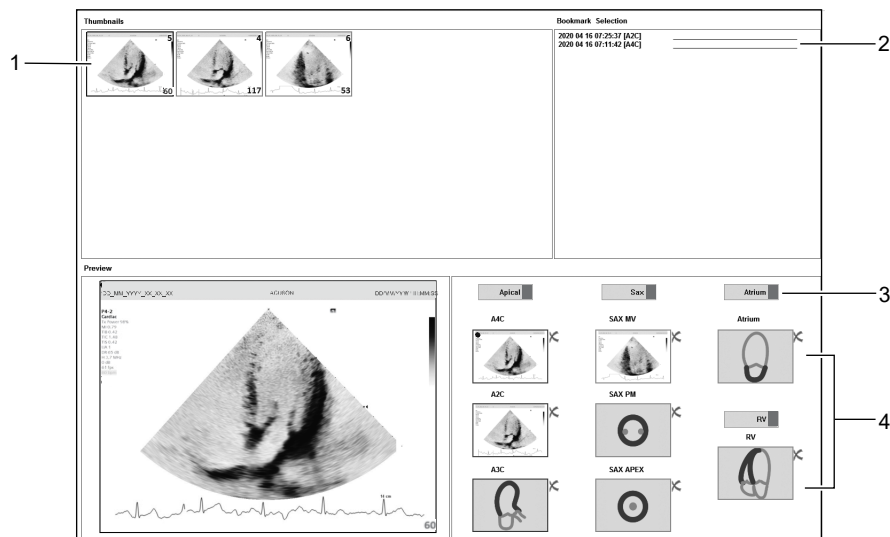

Example of the syngo VVI screen with multiple clips selected.

- 1 Thumbnails of clips selected during review
- 2 Bookmarks for accessing previously analyzed views
- 3 View selection labels
 - **Apical**
 - **SAX** (short axis)
 - **Atrium**
 - **RV** (right ventricle)
- 4 Thumbnails of clips within the selected view, for example, A4C

To analyze a view:

1. Activate review and then select a 2D-mode clip.
See also: For information about activating review, refer to Chapter 5 in the Instructions for Use.
2. To select multiple clips, select the required 2D-mode clips.
Note: You can select multiple images for analysis, but you can only process one view at a time.
3. Tap **syngo VVI**.
Note: You can activate the program only if you select a valid clip.
See also: For information about supported clips, refer to page D4-3.
4. Click an image in the upper left thumbnails box and then select the corresponding view (Apical, SAX, Atrium, or RV).
Note: You can select only one thumbnail for each view. If you add a new thumbnail, it replaces the existing thumbnail.
5. Repeat step 4 for each additional image for analysis.
Note: Each additional image must be within the same view.
6. To process a clip, select an image in the required view and then click the right arrow on the lower right of the screen.

7. To load a previously analyzed view from a bookmark:
 - a. Select the required bookmark from the list on the upper right of the *syngo* VVI screen.
 - b. Click the right arrow on the lower right of the screen.

8. To delete an image, click the **X** on the upper right of the image.
9. To exit the program, tap **Exit**. Or, press **2D**.

Analyzing a Different Clip in the Same View

Note: To analyze a different clip within the same view, you must have previously selected multiple clips for analysis. The selected clip displays as a thumbnail. You can select only one thumbnail for each view. If you add a new thumbnail, it replaces the existing thumbnail.

To process a trace for a different clip in the same view:

1. Select a clip in the active view.
2. Click the right arrow on the lower right of the screen.

Analyzing a Clip in a Different View

Note: To analyze a clip in a different view, you must exit *syngo* VVI, select an image in the required view, and then restart the program.

To begin processing a trace from a different view:


1. To exit the program, tap **Exit**.
2. Select a 2D-mode clip.
3. To select multiple clips from the review page, select the required 2D-mode clips.
4. Tap **syngo VVI**.
5. Click a clip in the upper left thumbnails box and then select the corresponding view (Apical, SAX, Atrium, or RV).
6. To begin processing a clip, select a clip in the required view, and then click the right arrow on the lower right of the screen.

Setting Up a Clip for Processing

If you load a clip with an ECG trace and two R-waves per heart cycle, *syngo* VVI defines the trace and performs the analysis.

If you load a clip without an ECG trace or without at least two R-waves per heart cycle, you must edit the clip margins, set R-wave markers, and adjust the end diastole and end systole frames.

Note: When a clip is loaded for analysis, ensure the end diastole and systole markers on the setup screen align with the ECG trace.

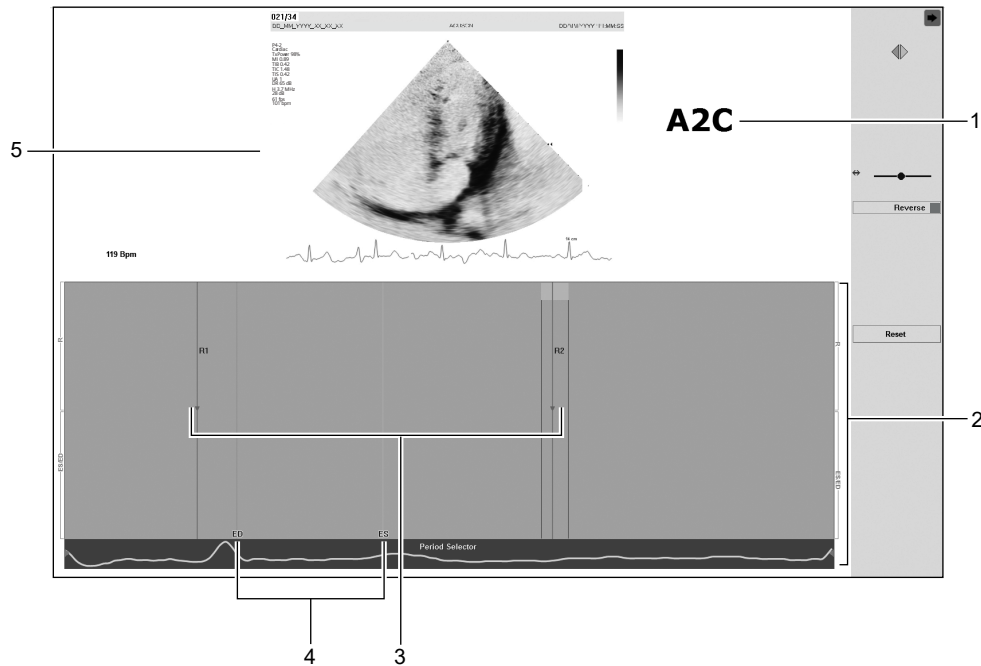

Example of the Setup screen.

- 1 Selected view of the heart
- 2 Period Selector and display
- 3 Current frame marker and R-wave markers
- 4 End diastole and end systole markers
- 5 2D clip with an ECG trace

Displaying M-mode Data

You can define an M-mode background for display behind each curve plot. The user-created drawing for the M-mode background can include multiple segments. An M-mode background provides additional timing information.

To display M-mode data by drawing an M-mode trace:

1. Click the sequence/M-mode selection button.



2. To display the M-mode sweep by drawing a trace:

- a. Position the cursor on the 2D image.

The system displays the pencil icon.



- b. Press **SET** at the beginning of the trace and then roll the trackball to create the trace.
 - c. To end the trace, press **SET**. Or, press **UPDATE**.

Note: It is not always necessary to draw lines. You can also use the CINE controls or trackball keys on the control panel to locate ECG R-waves and place markers as required.

3. Use the CINE controls to adjust the M-mode marker position corresponding to 2D-mode.

4. To add an R-wave marker:

- Click the required location in the upper half of the M-mode display.
- Double-click the required location in the lower half of the M-mode display and then click **Set Rwave**.

Note: ED-ES (end diastole-end systole) lines display in the lower half of the M line when you identify one complete cycle through the placement of R-waves.

5. To delete an R-wave marker:

- a. Position the cursor over a marker.
The system displays a double-line arrow.
 - b. Double-click the R-wave marker. Or, press **UPDATE**.

6. To reposition the markers for valve timing events:

- a. Position the cursor over a marker.
The system displays a double-line arrow.
 - b. Click the required marker and then roll the trackball to drag the marker to the required location.
 - c. Press **SET** to confirm the new location.

7. To reverse the grayscale of the M-mode sweep, for example, from black-on-white to white-on-black, select the **Reverse** check box.

8. To reverse all changes, click **Reset**.

9. To display the previous window, click the right arrow on the upper right of the screen.



Drawing and Modifying a Contour

The program displays the Contour screen if you load a clip with an ECG trace and R-wave information. Use the Contour screen to draw or edit an outline of the structure for analysis. You can display the Contour screen to modify or redraw the contour.



Example of the Contour screen.

- 1 Selections for displaying and editing the contour
- 2 Gamma slider
- 3 Placement guide for positioning the trace markers
- 4 Controls for repositioning the trace markers
- 5 ECG trace processed by syngo VVI (green line)
- 6 ECG trace displayed in the captured clip (gray line)
- 7 First marker indicated by a red dot
- 8 Heart rate
- 9 Imaging parameters and settings for the displayed clip transferred from the echo image

To create a trace:

Note: You can select a different image from the view on the lower right of the image screen.

1. Position the cursor over the required anatomy on the image.
The cursor displays as a +.

2. To create a trace using the placement guide:

See also: For an example of the placement guide, refer to page D4-9.

- a. Press **SET** to anchor the first marker.

A red dot displays at the selected point.

Target	Anchor location of first marker
Left ventricle	Base of the interventricular septum and mitral annulus
Left atrium	Mitral annulus near the interatrial septum
Right ventricle	Tricuspid valve annulus near the lateral wall

- b. Follow the placement guide on the screen to place subsequent markers.
 - c. Press **SET** to anchor each marker.
 - d. To anchor the final marker and complete the trace, press **SET** twice. Or, press **UPDATE**.
3. To manually create a trace:
 - a. Position the cursor over the next position for each segment of the trace, following a clockwise pattern.
 - b. Press **SET** to anchor each marker.

Note: Anchor 7 to 11 markers for apical traces. Anchor at least one marker to each wall segment. For short axis traces, do not anchor the final marker at or beyond the position of the first marker.
 - c. To anchor the final marker and complete the trace, press **SET** twice. Or, press **UPDATE**.
4. To reposition a marker, adjust the large markers and then adjust the subsequent smaller markers.
5. To change the brightness and contrast of the image:
 - a. Drag the Gamma slider to the right or left to adjust the gamma image setting.



- b. To restore the original brightness and contrast setting, click the restore gamma setting button.



6. To reverse the horizontal orientation of the image, click the pictures mirroring button.



7. To create a single endocardial trace, clear the **Endo+Epi** check box.
8. To create a simultaneous endocardial and epicardial trace:
- Select the **Endo+Epi** check box.
 - Click the green arrows on the right of the image to increase or decrease the distance between the endocardial and epicardial traces.




9. To process the trace, choose a method:
- Double-click the trace.
 - Press **UPDATE**.
 - Click the start analysis button.



A message indicates the analysis process is in progress.

Modifying a Contour

The system displays the contour for review and modification of the trace.

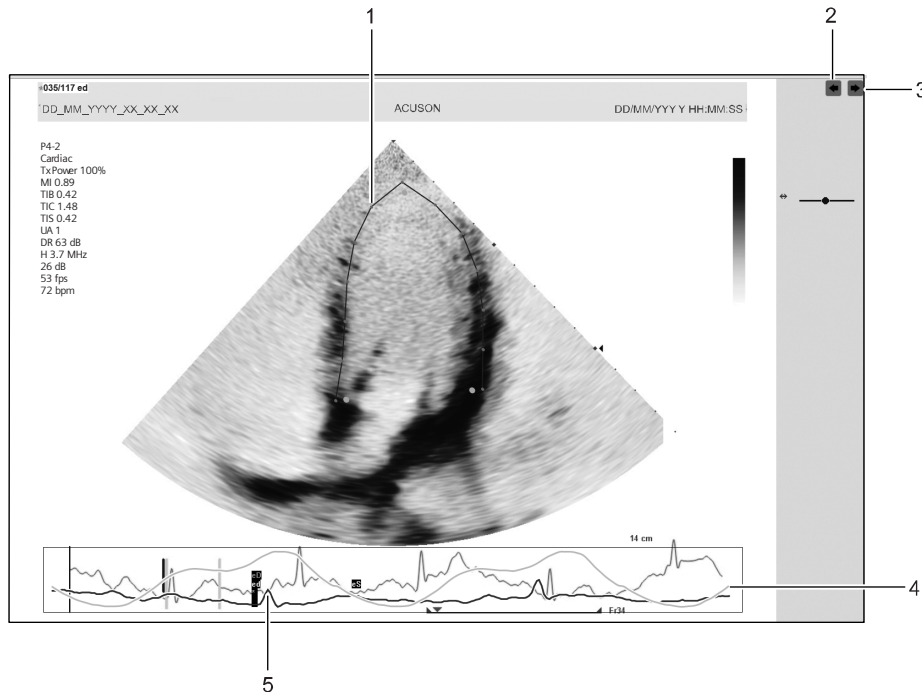

Example of a contour.

- 1 Drag the marker to reposition the marker and modify the contour.
- 2 Click the arrow to reject (delete) changes to the traced contour.
- 3 Click the arrow to accept changes to the traced contour.
- 4 Time volume curve
- 5 ECG trace

To modify a contour:

1. Use the CINE controls to review the trace.
2. To reposition a marker, adjust the large markers and then adjust the subsequent smaller markers.
3. To accept the changes to the contour and view analysis results, choose a method:
 - Double-click the trace.
 - Press **UPDATE**.
 - Click the accept changes arrow.



4. To reject the changes to the contour, click the reject changes arrow.




Processing a Contour and Viewing Results

When you process a contour, the program illustrates tissue velocity using velocity vectors (arrows). For each clip frame, the vector length illustrates the magnitude of velocity (speed), and the vector arrowheads illustrate the direction of motion.


Quantitative data derived from the velocity vectors are depicted as curve plots, parametric M-mode graphs, and parametric segmental diagrams. The specific results displayed depend on the view used to draw the contour, for example, apical view, short axis view, atrium view, or right ventricle view, and the selected parameter, for example, velocity, displacement, strain, or strain rate.

See also: For information about the selections for viewing analysis results, refer to page D4-18.

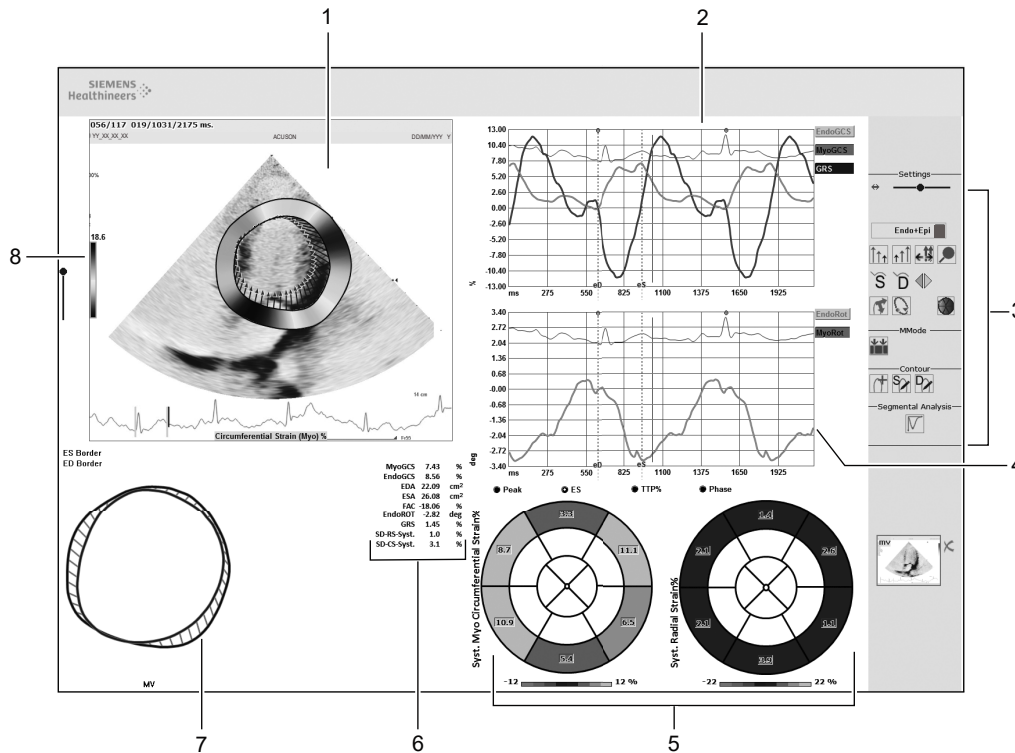

Example of the Analysis screen with clips displayed in SAX view.

- 1 2D clip with contour
 - 2 Global circumferential strain in the endocardial layer
 - 3 Selections for viewing analysis results
 - 4 Rotation of endocardial layer
 - 5 Systolic endocardial circumferential strain
 - 6 Analysis results for the selected clip
 - 7 End-systolic border and end-diastolic border
 - 8 Color bar for the selected parameter
- The system assigns a unique color bar to each parameter.

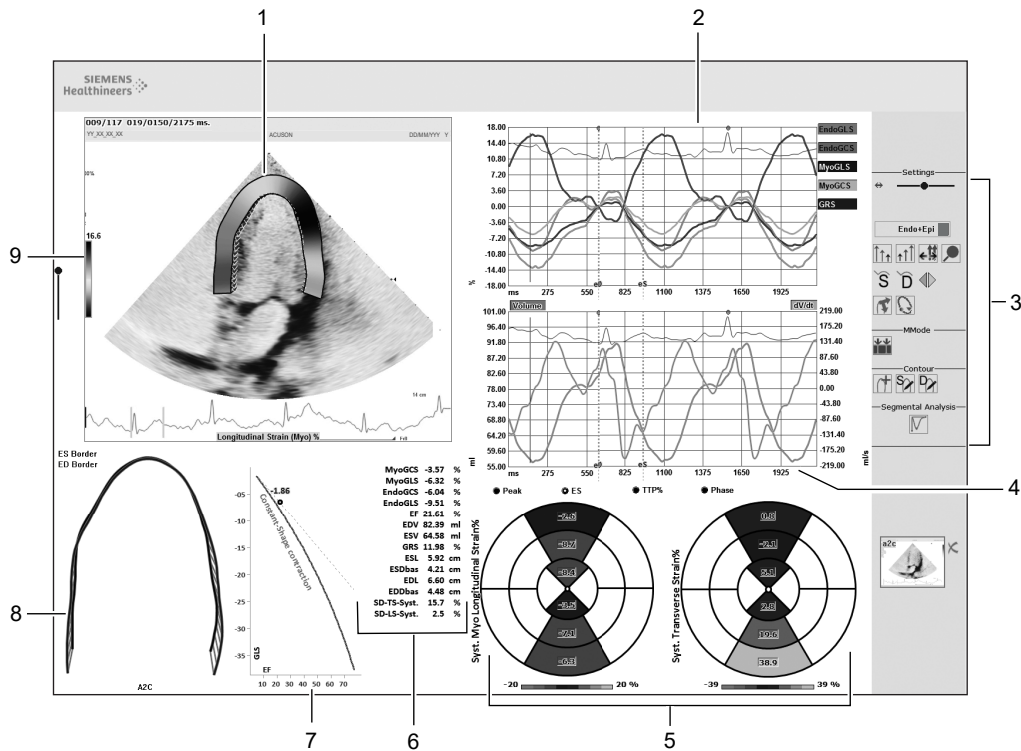

Example of the Analysis screen with clips displayed in Apical view.

- 1 2D clip with contour
 - 2 Global longitudinal and circumferential strain in the endocardial layer
 - 3 Selections for viewing analysis results
 - 4 Rotation of endocardial layer
 - 5 Systolic endocardial longitudinal strain
 - 6 Analysis results for the selected clip
 - 7 Graph displaying the relationship between the left ventricle ejection fraction and left ventricle global longitudinal strain
 - 8 End-systolic border and end-diastolic border
 - 9 Color bar for the selected parameter
- The system assigns a unique color bar to each parameter.

Viewing Results in a Bullseye Plot

(Available only for left ventricle images, not available for right ventricle or atrium images)

The bullseye plots support a 16-segment model. You can display summarized results of all analyzed views in the following segmental display formats (bullseye plots):


- **Overall Peak or Systolic Peak**

Note: The program displays overall peak results for the strain and displacement parameters and the systolic peak for the strain rate and velocity parameters.

- **Time to Peak**

The bullseye plots contain results from all clips analyzed in the current analysis group (apical or short axis). The Analysis step displays results only for a particular image. When you perform an analysis for one view of a patient study, the program adds the information to the bullseye plots. When you review an existing study, you must load at least one previously analyzed clip to update the bullseye plot.

Note: If you deactivate a segment in the Analysis step, the results for that segment are not included in the bullseye results and overall average.

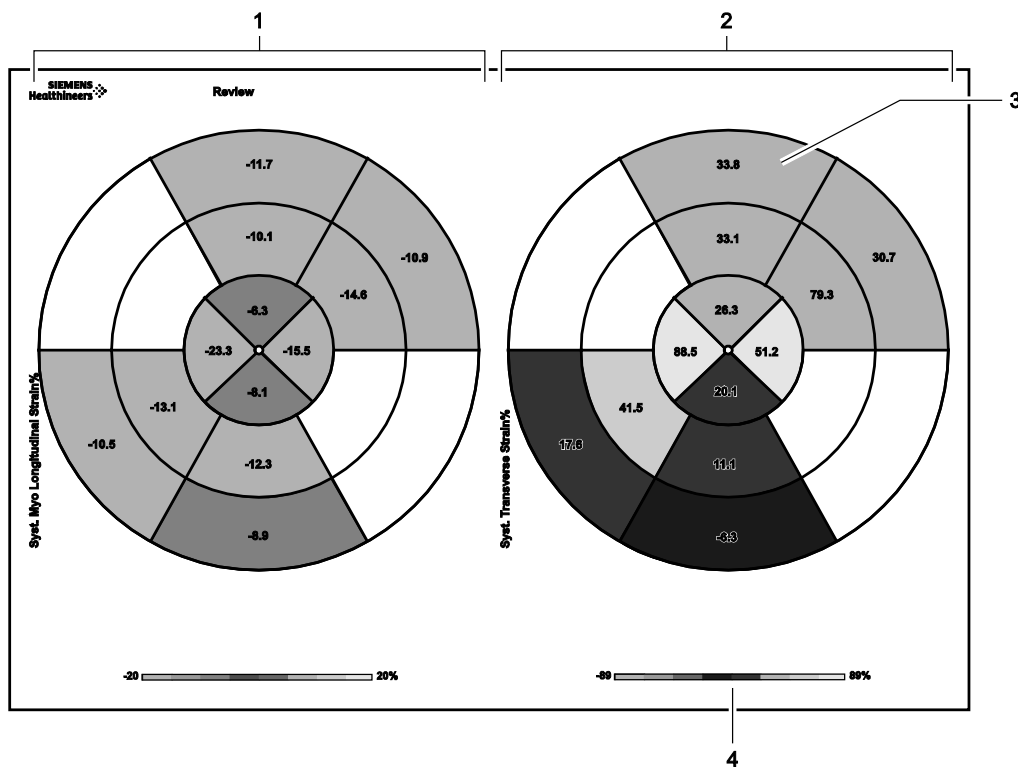

Example of bullseye plots and data.

- 1 Bullseye plot
- 2 Bullseye plot for strain results
(Available only when the **Endo+Epi** selection is enabled.)
- 3 Numerical result for the segment of the selected parameter
- 4 Color bar for the selected parameter; a unique color bar is assigned to each parameter

Analysis Page: Contour

Editing and Deleting a Trace

You can edit systole and diastole frames for a processed trace (contour) after you have completed creating the trace. After editing the trace on end systolic frames for the left ventricle or right ventricle and end diastolic frames for the left atrium, perform the analysis again.

To edit the currently displayed trace:

1. To modify the trace, click the required button.
 - To delete the current trace and starts a new trace, click the new trace button.

 - To modify both end systole and end diastole borders, click the end systole correction button.

 - To modify the end diastole border only, click the end diastole correction button.

2. Create a new trace or modify the displayed trace.
See also: For information about creating or modifying a trace, refer to page D4-9.
3. To process the modified trace, choose a method:
 - Double-click the trace.
 - Press **UPDATE**.
 - Click the start analysis arrow.

4. To view analysis results, choose a method:
 - Double-click the trace.
 - Press **UPDATE**.
 - Click the accept changes arrow.


To add endocardial and epicardial borders:

1. Select the **Endo+Epi** check box.
2. Click the green arrows on the right of the image to increase or decrease the distance between the endocardial and epicardial traces.



3. To process the modified trace, choose a method:
 - Double-click the trace.
 - Press **UPDATE**.
 - Click the start analysis button.



4. To view analysis results, choose a method:
 - Double-click the trace.
 - Press **UPDATE**.
 - Click the accept changes arrow.



Analysis Page: Settings

To	Do This
Change the brightness and contrast of the image	- Drag the gamma slider to the right or left to adjust the gamma image setting.  - To restore the original brightness and contrast setting, click the restore gamma setting button. 
Adjust the distance between the endocardial and epicardial traces	- Select the Endo+Epi check box. - Click the green arrows on the right of the image to increase or decrease the distance between the endocardial and epicardial traces. 
Activate or deactivate the zoom function	Note: When you use the zoom function, the system displays the image in a 2D full-screen. Click the zoom in/out button. 
Advance to the end systole	Click the go to end systole button. 
Advance to the end diastole	Click the go to end diastole button. 
Reverse the orientation of the image	Click the pictures mirroring button. 
Invert color parametric	Click the switch parametric over B-mode button. 
Set the position between the anterior and anteroseptal segment	(Available only in SAX view) Click the position between anterior and anteroseptal segment button. 

To	Do This
Display controls to toggle the contour, vector, orbit line and 2D of the trace	Click the toggle button. 
Change the size of the vector	Prerequisite: Display vectors before using these controls. - To decrease the height of the vector, click the decrease vector size button.  - To increase the height of the vector, click the increase vector size button.  - To reset the vector to its original size, click the reset vector size button. 

Analysis Page: Display Buttons

Note: The display buttons are available only when the toggle button is selected.

To view the display controls:

- Click the toggle button.



To	Do This
Display a color overlay with vector maps for regional strain	Click the arrows with color button.
Display vector maps for regional strain	Click the arrows button.
Display a color overlay with the contour	Click the borders with color button.
Display the contour only	Click the borders button.
Display orbit lines	Click the trajectories button.
Display segmentation deformation	(Available only when the epicardial border is displayed) Click the segmental deformation button.
Hide all overlays	Click the no overlay button.

Analysis Page

M-Mode

To	Do This
Select M-mode	Click the sequence/M-mode selection button. 

Contour


To	Do This
Begin a new trace	Click the new trace button. 
Correct the end systole border	Click the end systolic border correction button. 
Correct the end diastole border	Click the end diastolic border correction button. 

Segmental Analysis

To	Do This
Begin segmental analysis	Click the segmental analysis button. 

Segmental Analysis Results

The segmental analysis page displays peak and timing information, for example, time-to-peak values related to velocity, displacement, strain, and strain rate.

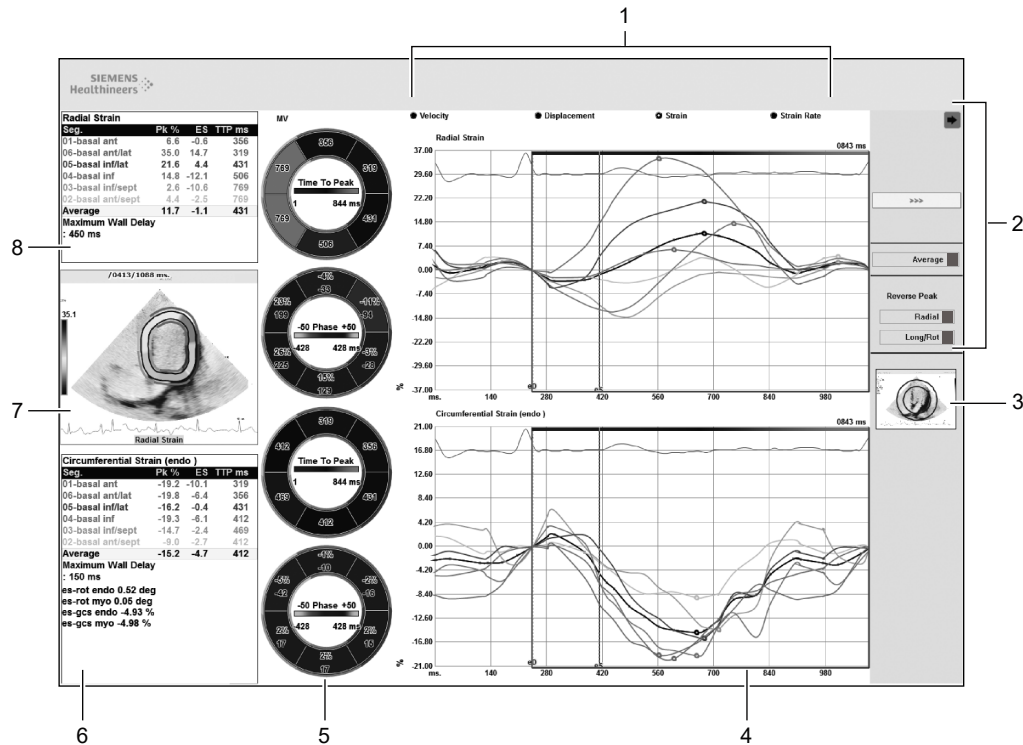

Example of the Segmental Analysis Page.

- 1 Display controls for velocity, displacement, strain, and strain rate
Strain is selected in this example.
- 2 Display and editing selections
- 3 Layer selection area
- 4 Circumferential strain information
- 5 Time to peak and phase information
- 6 Summary of circumferential or longitudinal strain data
- 7 Contour
- 8 Summary of radial or transverse strain data

Viewing Segmental Analysis Results

Use segmental analysis to view segmental curves and analysis results for segments of the contour. The program displays the clip with the contour, segment diagrams with time-to-peak and phase results, and curve plots for the selected segments. For apical views, the segmental analysis includes a rocking value to indicate motion of the apex.

- To activate or deactivate segments on the segment model, click the required segment on the segment model.

The curve plots display results only for the activated segments on the segment model.

- To display results for a selected location on the curve plot, position the pointer on the required data point on the curve.

The program displays the numerical results, segment name and number, and time for the selected location.

- To display or hide the average curve on the curve plots, click **Average**.

The average curve depicts the average of the analysis results for the currently displayed curves.

- To display the radial results:

- For minimum radial results, select the **Reverse Peak - Radial** check box.
- For maximum radial results, clear the **Reverse Peak - Radial** check box.

- To display the longitudinal/rotational results:

- For minimum longitudinal/rotational results, select the **Reverse Peak - Long/Rot** check box.
- For maximum longitudinal/rotational results, clear the **Reverse Peak - Long/Rot** check box.

- To display analysis results by parameter, click the required display format:

- **Velocity**
- **Displacement**
- **Strain**
- **Strain Rate**

- To display a different layer of the heart wall, click the required layer:

Note: Available only when the **Endo+Epi** selection is enabled.

- **Endo** displays the endocardial layer.
- **Myo** displays the myocardial layer.
- **Epi** displays the epicardial layer.

Analysis Display Information

Control	Description
Reverse Peak	For all parameters except strain longitudinal, strain circumferential, strain rate longitudinal, and strain rate circumferential, identifies the lowest point, the point that is 75% of the lowest point, or the point that is 50% of the lowest point, as determined by the "Time To" selection. For parameters strain longitudinal, strain circumferential, strain rate longitudinal, and strain rate circumferential, identifies the highest point, the point that is 75% of the highest point, or the point that is 50% of the highest point, as determined by the "Time To" selection.
	Moves from beat to beat in a continuous cycle when clicked. Disabled if only one beat is captured.
Velocity	Computes and displays values for the selected parameter. ▪ Longitudinal Velocity in Apical, Atrial, and Right Ventricular views ▪ Rotation Velocity is the rotation rate in SAX view. Counterclockwise rotation is positive. ▪ Long/Rot is the longitudinal velocity/rotation velocity rate in SAX view. Counterclockwise rotation is positive. ▪ Radial Velocity: – Transverse Velocity in Apical, Atrial, and Right Ventricular views – Radial Velocity in SAX view
Displacement	Computes and displays values for the selected parameter. ▪ Longitudinal Displacement in Apical, Atrial, and Right Ventricular views ▪ Rotation Displacement in SAX view ▪ Long/Rot is the longitudinal velocity/rotation rate in Apical view. Counterclockwise rotation is positive in SAX view. ▪ Radial Displacement requires an endocardial and epicardial trace. – Transverse Displacement in Apical, Atrial, and Right Ventricular views

Control	Description
Strain	Computes and displays values for the selected parameter. ▪ Longitudinal Strain in long axis and generic curve Apical, Atrial, and Right Ventricular views ▪ Long/Rot is the longitudinal/rotation rate in SAX view. Counterclockwise rotation is positive. ▪ Circumferential Strain in SAX view ▪ Radial Strain requires an endocardial and epicardial trace. – Transverse Strain in Apical and Right Ventricular views – Radial Strain in SAX view
Strain Rate	Computes and displays values for the selected parameter. ▪ Longitudinal Strain Rate in Apical, Atrial, and Right Ventricular views ▪ Long/Rot is the longitudinal velocity/rotation rate in SAX view. Counterclockwise rotation is positive. ▪ Circumferential Strain Rate in SAX view ▪ Radial Strain Rate requires an endocardial and epicardial trace. – Transverse Strain Rate in Apical and Right Ventricular views – Radial Strain Rate in SAX view

Saving Analysis Results to a Report

You can export analysis results to the patient report.

See also: For information about reports, refer to Chapter C1 in this manual.

To save data to a report:

- Tap **Save to Report**.

Viewing Analysis Results

The specific results displayed on the Analysis screen and exported to the worksheets and report depend on the view used to draw the contour for Apical view or SAX view and the selected parameter, for example, velocity, displacement, strain, or strain rate.

Results by Parameter

The Summary screen and worksheets display radial and longitudinal or circumferential results for each heart wall segment for the selected parameter: velocity, displacement, strain, and strain rate.

The Summary screen also displays the following summary results for the selected parameter:

- Average of all heart wall segments
- Standard deviation
- Maximum opposing wall delay in Apical view
- Maximum wall delay in SAX view

Exporting Images, Clips, and Data

You can export data in the following file formats to an external USB-compatible storage device or a network location.

- ASCII data in a TXT file
- XML data in a Microsoft Excel file

You can also export images and clips to an external USB-compatible storage device or a network location. Exported image files contain the displayed window without the on-screen selections. Exported clip files contain the clip in motion, as displayed within the clip window.

To export data or images:

1. To export data, tap **Export Data**.
2. To export images, select the required images and then tap **Export Image**.
3. Select the available media from the **Destination** list.
4. Click **Export**.

A message indicates the export is in progress.

See also: For descriptions of exported data files, refer to Appendix H in the System Reference.

To export a clip:

Note: You can export a clip as a video clip. You cannot export a clip as a DICOM clip.

1. Select the required clip.
2. To change the clip frame, use the CINE controls.
3. Tap **Export Clip**.
4. Select the available media from the **Destination** list.
5. Click **Export**.

A message indicates the export is in progress.

ACUSON Juniper
ACUSON Juniper Select
Diagnostic Ultrasound System
System Reference



An updated User and Reference Manual may be available for your ultrasound system. Access the most recent manual from the Siemens Healthineers Document Library at doclib.siemens-healthineers.com.

Product Version 3.0

Software Version VC10

©2018-2025 Siemens Medical Solutions USA, Inc.

All Rights Reserved.

Date of first issue: 2025-08

All trademarks are the property of their respective owners.

Siemens Medical Solutions USA, Inc. reserves the right to change its products and services at any time. In addition, this publication is subject to change without notice.

Manufacturer

Siemens Medical Solutions USA, Inc.
Ultrasound
22010 S.E. 51st Street
Issaquah, WA 98029
U.S.A.
Phone: +1-888-826-9702
siemens-healthineers.com

Siemens Healthineers Headquarters

Siemens Healthineers AG
Siemensstr. 3
91301 Forchheim
Germany
Phone: +49 9191 18-0
siemens-healthineers.com

Contents

Chapter 1	General System Configuration Description of the selections for customizing the ultrasound system using the configuration menu. Includes a procedure for backing up and restoring configuration settings and user-defined exams.
Chapter 2	Network Configuration Description of the selections for configuring network settings for the ultrasound system and the devices for storing and retrieving patient information, including worklist servers and network servers.
Chapter 3	Measurement and Report Configuration Description of the selections for customizing measurements and calculations for all exam types using the configuration menu.
Chapter 4	Clinical Applications Configuration Description of the selections for customizing clinical applications, for example Stress Echo, using the configuration menu.
Chapter 5	Export and Import Procedures for exporting and importing studies, including teaching files and slide shows.
Chapter 6	Security Configuration Information for the ultrasound system administrator. Includes procedures for creating user accounts and configuring the security log file.

Appendix A	EMC Publication (Electromagnetic Emissions and Immunity: Guidance and Manufacturer's Declaration) Information regarding the electromagnetic compatibility (EMC) testing of this ultrasound system.
Appendix B	Measurements and Calculations Reference Description of general measurement and calculation methods, and results.
Appendix C	Obstetrical Reference Description of measurement and calculation labels for the obstetrical exam. Includes obstetric report reference information.
Appendix D	Cardiac Reference Description of measurement and calculation labels for the cardiac exam. Includes cardiac report reference information.
Appendix E	Exam-specific Reference Description of measurement and calculation labels for all exam types except obstetrical and cardiac exam. Includes exam-specific report information.
Appendix F	Obstetrical Bibliography Bibliography of references for obstetrical examinations.
Appendix G	Cardiac Bibliography Bibliography of references for cardiac examinations.
Appendix H	Data File Descriptions Description of data files for the following features: ▪ <i>syngo</i> Velocity Vector Imaging

Note: Not all features and options described in this publication are available to all users. Please check with your Siemens Healthineers representative to determine the current availability of features and options.

See also: For a list of available options and features, refer to Appendix A in the Instructions for Use.

About the User and Reference Manuals

The user and reference manuals contain descriptions for the following ultrasound systems:

- ACUSON Juniper diagnostic ultrasound system
- ACUSON Juniper Select diagnostic ultrasound system

Features and options unique to an ultrasound system are identified within the following locations in the user and reference manuals:

- Chapter 1, Instructions for Use
- Chapter 3, Instructions for Use
- Appendix A, Instructions for Use

The user and reference manuals consist of the following publications.

Publication	Includes
<i>Instructions for Use</i>	▪ Conventions and typographical conventions used in the manuals ▪ Intended audience ▪ Introduction to the ultrasound system ▪ Technical description of the ultrasound system ▪ Safety and care information for the system and compatible transducers ▪ Descriptions of system controls ▪ Procedures for system setup, examination fundamentals, image acquisition and optimization, and examination review ▪ Information about system security and protection of patient data against security threats ▪ Acoustic output data
<i>Clinical Applications Guide</i>	▪ Explanation of clinical applications and clinical software programs for use on the ultrasound system ▪ Procedures for the physiologic function and the biopsy function ▪ Description of the transesophageal transducer, including cleaning and care of the transducer ▪ General and exam-specific measurements and calculations ▪ Exam-specific annotations
<i>System Reference*</i>	▪ General and exam-specific measurement and calculation labels ▪ Description of customizable system settings ▪ Information about DICOM connectivity and network capabilities ▪ Clinical references ▪ Information regarding the electromagnetic compatibility (EMC) testing of this system

*Languages supported by the user interface include a translation of this publication.

Conventions

Take a moment to familiarize yourself with the following conventions in your user manual.

Convention	Description
Warnings and Cautions	 WARNING: Warnings are intended to alert you to the importance of following the correct operating procedures where risk of injury to the patient or system user exists.  Caution: Cautions are intended to alert you to the importance of following the correct operating procedures to prevent the risk of damage to the ultrasound system.
Notes and Prerequisites	Prerequisite: Prerequisites contain tasks the user must complete or information the user needs prior to performing a procedure. Note: Notes contain information concerning the proper use of the ultrasound system or correct execution of a procedure.
Examples of Cross-References	Examples: See also: For information about configuring a network, refer to Chapter 2 in the System Reference. See also: Touch screen controls, page 1-17 See also: For a list of transducers compatible with the needle guide accessories, refer to Chapter 6 in this manual.

Typographical Conventions

The following typographical conventions are in descriptions and procedures within your user and reference manuals. Use the conventions to identify the location of a control or selection on the ultrasound system.

Control or Selection	Description
Controls on the Control Panel and Keys on the Keyboard	Controls on the control panel are indicated by uppercase, boldface type. The function of the trackball and trackball keys is indicated by on-screen labels displayed on the image screen. Examples: ▪ Push DEPTH / FOCUS. ▪ Rotate 2D to adjust the 2D gain. ▪ Press CALIPER to activate the measurement function. ▪ Press and hold the Alt key on the keyboard to access the special characters. ▪ Press UPDATE. ▪ Position the marker and then press SET.
Controls on the Control Panel and Soft Key Selections on the Touch Screen	The touch screen has unlabeled controls located on the control panel. Tap, press, rotate or press and then rotate the control to make a selection or change a setting. The selection assigned to these controls is indicated by boldface type and enclosed in square brackets ([]). Examples: ▪ Press [Tx Power] downward to decrease transmit power. ▪ Rotate [Map].
Selections on the Image Screen	On-screen selections are indicated by boldface type. Roll the trackball to position the pointer on the on-screen control or object and then press a trackball key. Examples: ▪ Click Cancel. ▪ Double-click a thumbnail. ▪ Select the required check box to search for scheduled procedures.
Selections on the Touch Screen	On-screen selections are indicated by boldface type. To access the selection, tap the touch screen. To change the direction of an object displayed on the touch screen, drag the object on the touch screen. Examples: ▪ Tap Stress Echo. ▪ Drag the blue dot to indicate the direction of the arrow. ▪ Swipe the page to select the next page.

1 General System Configuration

System Configuration Overview.....	2
Categories of Settings for System Configuration.....	2
General 1.....	3
General 2.....	4
Display.....	5
Custom Keys.....	6
Soft Keys.....	9
Procedure.....	10
Assigning an Exam Preset to a Procedure Description.....	10
Exam.....	11
Exam Association.....	13
Assigning Measurement Packages and Annotations to an Exam.....	13
Measurement.....	13
Text Annotation.....	14
Body Marker.....	15
Stress Echo.....	15
Protocols.....	15
Importing and Exporting Workflow Protocols.....	18
Contrast.....	19
Clip Store.....	19
Printer.....	20
Configuring and Installing Printers.....	20
Configuring Printing Preferences.....	21
Barcode.....	22
Wired and Wireless Selections.....	23
Configuring Wired Connections to a Network.....	23
Configuring Wireless Connections to a Network.....	23
Network.....	24
Utility.....	24
Archive.....	24
Security.....	25
Enabling Ports.....	26
Importing or Exporting a Client Certificate.....	26
Service.....	27
Enabling Anti-Virus Protection.....	27
Backing Up and Restoring Configuration Settings.....	27

System Configuration Overview

You can modify the configuration settings at any time.

Categories of Settings for System Configuration

Menu Item	Description
General 1	Configure system-wide settings, including the name of the hospital or clinic, system language, date and time formats, audio settings, trackball behavior, and shutdown preferences.
General 2	Configure settings for patient data, annotations, storing images, and S-video.
Display	Configure image screen settings, including screen saver behavior, display of the depth gain compensation, zoom navigation, and Doppler and M-mode cursors.
Custom Keys	Assign functionality to user-definable controls on the control panel and the footswitch.
Soft Keys	Customize the functionality assigned to soft keys and controls.
Procedure	Configure settings for defining a requested or scheduled procedure.
Exam	Configure settings that can vary by exam type, including freeze behavior, dual display format, physio management, and settings specific to M-mode, Doppler, or color imaging.
Exam Association	Configure settings for assigning a default measurement package, a library of body markers, and a library of text annotations to an exam.
Measurement	Select measurement conventions and establish measurement and report preferences for each type of exam.
Text Annotation	Customize libraries of text for each type of exam.
Body Marker	Customize libraries of body markers for each type of exam.
Stress Echo	(Available only for the cardiac exam type) Configure settings for Stress Echo workflow protocols.
Protocols	Configure settings for user-defined workflow protocols.
Contrast	Configure settings for capturing clips during contrast agent imaging.
Clip Store	Configure settings for capturing clips, including capture length, capture type, and compression.
Printer	Configure printers connected to the system.
Barcode	Configure settings to use a barcode reader to enter patient registration information.
Wireless	Configure settings for connection to a wired or wireless network.
Network	Configure network connections for the local host, storage server, worklist server, and network export server.
Utility	Configure settings for clip playback and slide shows.
Archive	Configure settings for the ultrasound system after studies are archived, including deletion of study history.
Security	- Indicates enabled or disabled security features on the ultrasound system. - Provides selections for configuring security settings.
Service	Provide access for authorized Siemens Healthineers service representatives to troubleshoot and maintain the ultrasound system.
Close	Save changes and exit the current page.

To modify the system configuration settings:

1. Tap the configuration icon.



The system displays the configuration menu on the left side of the screen and the settings for the selected menu option on the right.

2. Click the menu option for the required settings.
3. Enter the required settings.
4. Repeat steps 2 and 3 for each category of required settings.
5. To exit the configuration settings, click **Close**.

General 1

Selection	Description
Organization	---
Institution Name	Name of your hospital or clinic
Regional Settings	---
Language	System language ▪ English ▪ German ▪ French ▪ Spanish ▪ Italian ▪ Chinese ▪ Russian (applies only to the keyboard)
Date Format	Format for the system date: ▪ Month/Day/Year ▪ Day/Month/Year ▪ Year/Month/Day
Date	Current date
Time	Current time
Time Zone	Current time zone
Height and Weight Units	Units for the patient's height and weight: metric or U.S.
Audio	---
Beep on Button Press	Enable or disable an audible tone when a control is pressed
Beep after Store Completes	Enable or disable an audible tone when storage is completed
Beep Volume	Volume of the beep
Doppler Volume	Volume of Doppler audio

Selection	Description
Power Off	Response of the power on/off control: ▪ Shut Down with Confirmation Dialog displays a confirmation message when powering off the system ▪ Shut Down immediately powers off the system
Trackball Behavior	---
Trackball Speed	Select the responsiveness of the system to trackball movement. ▪ Low repositions an object a short distance with minimal trackball movement. ▪ Medium repositions an object a moderate distance with minimal trackball movement. ▪ Fast repositions an object a long distance with minimal trackball movement.
System Start Up	---
Active Transducer	Select a transducer port displayed as selected when you power on the ultrasound system. ▪ T1 ▪ T2 ▪ T3 ▪ T4 ▪ T5

General 2

Selection	Description
Patient	---
Store Patient Screen on Registration	Enable or disable storing a screen representation of the completed patient data form to the registered patient's study.
Hide Patient Banner	Show or hide the patient banner on the image screen.
Patient Data	Show the patient data either by age or by date of birth.
End-Exam Behavior	---
Activate the 2D imaging mode	Activate 2D imaging when the exam ends.
Register a patient or search from Worklist	Display the patient registration form to register a new patient or display the worklist to search for a patient study if the worklist server is configured.
Physician or Operator	Add or display the initials or other identifying information of the person performing and referring the patient study. ▪ Add registers the user-defined physician or operator name. ▪ Delete deletes the selected physician or operator name. ▪ (list) displays the user-defined physician or operator names entered in the configuration settings or patient registration form.

Selection	Description
Annotation	---
Delete Body Marker on Unfreeze	Show or hide the body markers when an image is unfrozen.
Delete Text on Unfreeze	Show or hide text and arrow annotations when an image is unfrozen.
Arrow Type	Select the type of the directional arrows for annotations.
Arrow Size	Select the size of the directional arrows for annotations.
Arrow Color	Select the color of the directional arrows for annotations.
Image Store	---
Store Image with Caliper	Enable or disable the display of measurement markers with a saved image.
S-Video	---
S-Video Format	Select the signal standard format for separate video connection, for example, National Television Standards Committee (NTSC) or Phase Alternating Line (PAL).

Display

Selection	Description
Monitor	---
Enable Screen Saver	Enable or disable the screen saver.
Screen Saver Delay	Time delay in minutes before the system activates the screen saver
Screen Saver Type	Select the background for the screen saver.
DGC	---
DGC Invert with Image Invert	Invert the graphic indicating the DGC curve with the image
DGC Curve Display	Select the display of the depth gain compensation curve on the image screen ▪ Off does not display the depth gain compensation curve. ▪ Always On displays the depth gain compensation curve continuously. ▪ Time Out temporarily displays the depth gain compensation curve after edits to the depth gain compensation.
DGC Control	Select the image depth setting, for example, the current image depth or maximum image depth
Zoom Navigation	---
Display Zoom Navigation on Image	Enable or disable the display of the zoom navigation for high-definition zoom on the image screen
Zoom Navigation Size	Specify the size of the zoom navigation for high-definition zoom ▪ Small ▪ Medium ▪ Large

Selection	Description
Doppler / M-Mode	---
Bypass Cursor Mode	Select the response of the M-mode or Doppler control: ▪ Immediately display the M-mode sweep or Doppler spectrum ▪ Display the M-mode or Doppler cursor in the 2D image and then requires a second press of the mode control to display the sweep or spectrum
Doppler Frequency/Velocity	Unit in the Doppler spectrum: frequency (Hz) or velocity (cm/s)
Doppler Auto Optimization	Select the parameters to use with Doppler auto optimization technology. ▪ Gain ▪ Dynamic Range ▪ Scale / Baseline ▪ Auto Invert

Custom Keys

Selection	Description
Key Direction	---
Increase Zoom	Assign a clockwise or counterclockwise rotation to the control to increase the magnification factor
Depth / Focus	Assign a clockwise or counterclockwise rotation or upward or downward movement to the control to adjust the imaging depth or focus: ▪ Move up Focus assigns a clockwise or counterclockwise rotation to adjust the imaging focus and an upward or downward movement to adjust the imaging depth. ▪ Increase Depth assigns a clockwise or counterclockwise rotation to adjust the imaging depth and an upward or downward movement to adjust the imaging focus.
Key Function	---
Print Store 1 Print Store 2	Assign the required functions to the selected control on the control panel: ▪ PC Printer ▪ Image Store ▪ Clip Store ▪ USB BW ▪ USB Color ▪ Image Store & DICOM Printer ▪ Image Store & PC Printer ▪ Image Store & USB BW ▪ Image Store & USB Color
Clip Store	Assign one of following functions to the clip store key. ▪ Image Store ▪ Clip Store

Selection	Description
Left Pedal	Assign one of the following functions to each pedal on the footswitch: ▪ Freeze ▪ Image Store ▪ Clip Store ▪ Caliper ▪ Update ▪ TEQ Update ▪ U/D Flip ▪ L/R Flip ▪ USB BW ▪ USB Color ▪ Rotate +90° (TEE) ▪ Rotate -90° (TEE) ▪ Center Line ▪ Live Bi-plane ▪ Toggle View
Right Pedal	

Selection	Description
User Defined 1	Assign one of the following functions to the selected control on the control panel:
User Defined 2	▪ Full D activates a full-screen Doppler spectrum
User Defined 3	▪ Full M activates a full-screen M-mode sweep ▪ Arrow activates the arrow annotation function ▪ Clip Store acquires and saves a clip or static image ▪ Biopsy activates the biopsy function ▪ TEQ Update applies updated TEQ optimization settings to the image ▪ << Steer adjusts the angle of the field of view to the left for linear array transducers ▪ Steer >> adjusts the angle of the field of view to the right for linear array transducers ▪ Text A activates the library of annotations as specified in the configuration settings ▪ Text B activates the library of annotations as specified in the configuration settings ▪ DVR enables or disables recording for the DVR connected to the system ▪ Trapezoid enables or disables a trapezoid image formats ▪ U/D Flip reverses the orientation of the image in vertical (up/down) ▪ L/R Flip reverses the orientation of the image in horizontal (left/right) ▪ Full Size resets the field of view to the full size of the image area ▪ Needle Visual. activates the needle visualization function ▪ Pointer activates the mouse cursor on the screen ▪ CW activates the continuous wave mode ▪ Simult enables or disables the simultaneous display of the 2D image and Doppler spectrum in real time ▪ Dual enables or disables the display of one real-time and one frozen image in a side-by-side format ▪ Live Dual enables or disables the display of two real-time images on the image screen ▪ Quad enables or disables the display of four separately acquired images ▪ None indicates no function is assigned to the control ▪ Delete deletes the measurement results, annotations, body markers, and the selected images or clips ▪ Rotate +90° (TEE) rotates the image for the transesophageal transducer clockwise in 90-degree increments ▪ Rotate -90° (TEE) rotates the image for the transesophageal transducer counterclockwise in 90-degree increments ▪ Center Line shows or hides the vertical dotted line depicting the center line of the image ▪ Live Bi-plane activates or deactivates biplane imaging during real-time imaging ▪ Toggle View alternates between the transverse and sagittal views during biplane imaging ▪ Exit Dual exits dual display format ▪ pSWE activates or deactivates Point Shear Wave Elastography (pSWE) ▪ UDFF activates or deactivates Ultrasound-Derived Fat Fraction (UDFF)
Configure END Key	Assign one of the following functions to the selected control on the retractable keyboard: ▪ Not Used disables the key on the keyboard. ▪ Close Study closes the current study.

Soft Keys

You can customize the touch screen controls available during a selected operating mode, including live and frozen imaging for each mode.

To customize the display of soft key selections:

1. Access the configuration settings and then tap **Soft Keys**.
2. Select the required label.
 - **2D Lv** specifies 2D mode with a real-time image
 - **2D Frz** specifies 2D mode with a frozen image
 - **C Lv** specifies color mode with a real-time image
 - **C Frz** specifies color mode with a frozen image
 - **P Lv** specifies power mode with a real-time image
 - **P Frz** specifies power mode with a frozen image
 - **M Lv** specifies M-mode with a real-time image
 - **M Frz** specifies M-mode with a frozen image
 - **D Lv** specifies Doppler with a real-time image
 - **D Frz** specifies Doppler with a frozen image
3. To assign a new soft key selection to a button:
 - a. Select the required button representing the soft key selection.
 - b. Select the required label from the drop-down list.
4. To remove a soft key selection:
 - a. Select the required button.
 - b. Select **Blank** from the drop-down list.

Procedure

When you assign exam presets to your procedure descriptions, the patient registration form displays the exam preset when you register a patient from a worklist server.

Assigning an Exam Preset to a Procedure Description

You can assign an exam preset to a requested or scheduled procedure.

You can export procedure descriptions for importing onto an equivalent ultrasound system or to edit the procedure descriptions on a personal computer.

To assign an exam to a procedure description:

Prerequisite: Obtain procedure descriptions from your PACS administrator.

1. Access the configuration settings and then tap **Procedure**.
2. Select the type of procedure:
 - **Requested Procedure**
 - **Scheduled Procedure Step**

To	Do This
Add a new procedure	▪ Enter a procedure description and then click Add. The system adds the procedure to the Procedure Description list.
Assign an exam preset to a procedure	▪ Select the required procedure and then select an option from the Exam Presets list. Note: For user-defined exam presets, the name of the base exam preset displays in the Based on list. See also: For information about creating a user-defined exam, refer to Chapter 4 in the Instructions for Use.
Remove a procedure description	▪ Select the required procedure description and then click Delete.
Export the procedures	1. Connect a USB-compatible storage device to the ultrasound system. 2. Click Export. 3. To update the list of available media, click Refresh. 4. Select the name of the USB-compatible storage device. 5. Edit the name of the file, if necessary, and then click OK. 6. Remove the USB-compatible storage device from the ultrasound system.
Import procedures	1. Connect the USB-compatible storage device with the procedures file to the ultrasound system. 2. Click Import. 3. To update the list of available media, click Refresh. 4. Select the name of the USB-compatible storage device. 5. Select the required file and then click OK.

Exam

Selection	Description
Exam	Configuration settings are available for each selected exam type
(exam type)	Exam type for customization
Freeze Behavior	Assign one of the following functions to the selected control on the control panel: ▪ Caliper activates the measurement function ▪ Cine activates the CINE function ▪ Text activates the text annotation function ▪ Arrow activates the arrow function ▪ Body Marker activates the body marker function
Seamless Dual	Enable or disable seamless display of dual
Dual Behavior on Freeze	Assign one of the following functions to the Dual mode when the image is frozen: ▪ Left Frozen-Right Live displays a real-time 2D image on the right and a frozen image on the left. ▪ ED-ES Split displays an end-diastolic frame and end-systolic frame simultaneously.
Dual Display Format	Specify the dual display format. ▪ Side by Side displays either real-time images or frozen images in a side-by-side format. ▪ Top-Bottom displays either real-time images or frozen images in a vertical format.
Auto TEQ	When enabled, activates Auto Tissue Equalization (Auto TEQ) optimization
2D/M & 2D/Doppler Display Format	Image screen layout when two imaging modes are active: ▪ 1/2 2D – 1/2 Trace uses one-half of the image screen to present the 2D image above the M-mode sweep or Doppler spectrum. ▪ 1/3 2D – 2/3 Trace uses one-third of the image screen to present the 2D image above the M-mode sweep or Doppler spectrum. ▪ 2/3 2D – 1/3 Trace uses two-thirds of the image screen to present the 2D image above the M-mode sweep or Doppler spectrum. ▪ Side by Side presents the 2D image side-by-side the M-mode sweep or the Doppler spectrum.
Doppler Update Style	Display of the 2D image during Doppler (frozen or live), or the interval for updating the 2D image: ▪ 2D-Frz / D-Lv displays 2D frozen and Doppler live ▪ 2D-Lv / D-Lv displays 2D live and Doppler live ▪ D-Lv / 2D-1s displays Doppler live and 2D update in one second intervals ▪ D-Lv / 2D-2s displays Doppler live and 2D update in two second intervals ▪ D-Lv / 2D-4s displays Doppler live and 2D update in four second intervals ▪ D-Lv / 2D-8s displays Doppler live and 2D update in eight second intervals
Doppler Auto Trace	Enable or disable the display of measurements and graphical traces to track mean and/or maximum velocities or frequencies
Doppler Trace Orientation	▪ Above Baseline displays a trace above the baseline ▪ Below Baseline displays a trace below the baseline ▪ Both displays a trace above and below the baseline

Selection	Description
Auto Invert of Color While Steering	Enable or disable inverting the color velocity scale when the region of interest is steered
Auto Invert of Color and Spectrum	Enable or disable inverting the Doppler spectrum and the color bar together in 2D-mode with Color / Doppler when the invert control on the control panel is pressed
2D Steer with Cursor	Enable or disable steering of the 2D image when the Doppler cursor or color region of interest is steered
Auto 2D FOV with Color Activation	Adjust the 2D-mode field of view (FOV) to the color region of interest (ROI).
Physio Management	The lead type for use with the physiologic function
3 Lead ECG	Enable or disable the display of the related trace when the selected exam is activated. The ECG trace is from the ECG input
Aux ECG	Enable or disable the display of the external ECG trace from the auxiliary input for the selected exam
Aux Signal	Enable or disable the display of the external signal from the auxiliary input for the selected exam
Respiration	Enable or disable the display of the respiratory trace for the selected exam
Indicate R-Wave	Enable or disable R-wave triggers during image acquisition or CINE and indicate poor detection of an R-wave
Lead Type I	Lead type convention for the right arm and the left arm
Lead Type II	Lead type convention for the right arm and the left leg
3D/4D Default Selection	Assign 3D imaging or 4D imaging to the selected control on the control panel
Bi-plane Default View	Specifies the default plane in a single display format during biplane imaging ▪ Transverse ▪ Sagittal

Exam Association

When you assign a measurement package, a library of text annotations, and a library of body markers to the exam, the exam activates with the configured settings.

Assigning Measurement Packages and Annotations to an Exam

You can assign a default measurement package, a library of body markers, and a library of text annotations to an exam.

To assign preferences to an exam:

1. Access the configuration settings and then tap **Exam Association**.
Note: The list includes user-defined exams added during imaging.
2. Configure your preferences for the required exam:
 - **Exam** displays the list of exams.
 - **Measurement** specifies the default measurement package for the exam.
 - **Text** specifies the default library of text annotations for the exam.
 - **Body Marker** specifies the default library of body markers for the exam.

Measurement

You can configure the general and exam-specific settings for the measurement function.

See also: For information on configuring the measurement function, refer to Chapter 3 in this manual.

Text Annotation

Each exam type has a library of text annotations and modifiers for anatomical structures, body positions, and acquisition planes.

Note: The system identifies modified text with an asterisk (*).

Selection	Description
Exam	Exam type for customization
Text A, B, C, D	Text for identifying body orientation, direction, or proximity for the selected exam
(text box)	Text box for adding or modifying user-defined modifier text for the selected exam
Available Text A, B, C, or D (tab)	Predefined text available for adding to the list of text annotations and modifiers for the selected exam.
Selected Text A, B, C, or D (tab)	List of predefined text annotations and modifiers for the selected exam.
Library	List of predefined text for the selected exam
(text box)	Text box for entering user-defined text or editing the selected predefined text
Available Annotations	Predefined text available for adding to the list of text annotations for the selected exam
Selected Annotations	Text annotations for the selected exam
Enable Auto Complete	Enable or disable the display of predefined text matching the first letters you enter on the image screen
Modify	Replace the selected text with the characters entered in the text box.
(right arrow)	Add the selected text to the list of selected text.
(left arrow)	Remove the selected text from the list of selected text.
(up arrow)	Reposition the selected text upward in the list.
(down arrow)	Reposition the selected text downward in the list.
Restore Defaults	Deletes the user-defined settings and restores the factory default settings

To customize text annotations:

- Access the configuration settings and then tap **Text Annotation**.

To	Do This
Insert text into the annotation library or a modifier library	- Enter new text in the text box. Or, select predefined text in the list of available text. - Click the right arrow.
Edit text in the annotation library or in a modifier library	- Select the text from the selected text field. - Edit the text in the text box. - Click Modify.
Remove text from the annotation library or a modifier library	- Select the text in the list of selected annotations or modifiers. - Click the left arrow.
Reset all annotations to the factory settings	Click Restore Defaults .
Change the sequence of annotations or modifiers	- Select the text in the list of selected annotations or modifiers. - To reposition the text upward, click the up arrow. To reposition the text downward, click the down arrow.

Body Marker

You can configure body markers for a specific exam to indicate the anatomy under evaluation and the orientation of the transducer.

Selection	Description
Exam	Exam type for customization
Available Body Markers	Body markers available to add to the specified exam.
Selected Body Markers	Body markers selected for the specified exam
(right arrow)	Add the selected body marker to the list of selected body markers.
(left arrow)	Remove the selected body marker from the list of selected body markers.
(up arrow)	Reposition the selected body marker upward in the list.
(down arrow)	Reposition the selected body marker downward in the list.
Restore Defaults	Deletes the user-defined settings and restores the factory default settings.

Stress Echo

You can configure stress echo settings for capturing and reviewing cardiac loops for a multi-phase stress echo protocol.

See also: For information on configuring stress echo settings, refer to Chapter 4 in this manual.

Protocols

A protocol is a predefined checklist to guide you through a clinical workflow.

Selection	Description
Protocols	Displays an expandable and collapsible list of protocols and associated views.
Add Protocol	Inserts a blank protocol at the bottom of the protocol list.
Add View	Inserts a blank view at the end of the list of views for the selected protocol.
Copy	Copies the selected protocol or view, assigns a name that includes the copied protocol or view name with a numeric extension, and adds the new protocol or view.
Delete	(Not available for factory default protocols) Removes the selected protocol or view.
Import	Copies a protocol from an external storage device to the ultrasound system.
Export	Copies a protocol from the ultrasound system to an external storage device.
Restore Defaults	(Not available with user-defined protocols) Restores the factory default settings of the selected protocol.
Protocol Name	Displays the name of the selected protocol.
Progression	Predefines the next view in the protocol: ▪ Automatic activates the next view in the protocol when the current view is complete. ▪ Manual waits for the user to select the next view in the protocol when the current view is complete.
Automatically ends when all views completed	Automatically ends the protocol when the last view is complete.

Selection	Description
End Protocol Behavior	---
Remain in Acquisition	Continues acquisition when the selected protocol is ended.
Enter Review	Activates review when the selected protocol is ended.
Enter Report	Displays the patient report when the view of the selected protocol is ended.
View (tab)	Displays the current settings for the selected view.
View Name	Displays or modifies the name of the selected view.
Mode Selection	Lists available imaging modes to activate for the selected view.
D	Configures the Doppler settings: ▪ PW specifies pulsed wave Doppler ▪ CW specifies continuous wave Doppler ▪ DTI specifies Doppler tissue imaging
Color	Configures the color mode settings: ▪ Color selects color mode ▪ Power selects power mode ▪ DTI specifies Doppler tissue imaging ▪ DTE specifies Doppler tissue energy
2DMode	Configures the 2D-mode settings: ▪ Dual specifies dual 2D-mode ▪ Live Dual displays two real-time images from the same acquisition on the image screen
M-Mode	Selects M-mode.
Measurements (tab)	Displays the current measurement settings for the selected view.
Label1, Label2, or Label3 (tab)	Selects the first measurement label.
Application	Lists the available exams.
Imaging Mode	Selects the imaging mode to activate for the selected measurement.
Package / Group	Selects the measurement group.
label type (option)	Selects the measurement label for the selected measurement. ▪ Select a measurement label from the list and an additional measurement prefix, for example, right, left or none. ▪ Select an exam type from the list and an additional measurement prefix, for example, proximal, middle, distal, or none.
Press Caliper to Resume	(Available only for Label2 and Label3) Enables or disables the measurement function when the user presses the caliper control on the control panel.

Selection	Description
Text & Body Marker (tab)	Displays the current body marker and annotation settings for the selected view.
Text Annotation	Displays the on-screen annotation.
Location	Provides settings for defining the default location of the text annotation. ▪  Previews the user-defined location on the image screen for text annotation ▪  Defines the default location on the image screen for the text annotation
Apply to all Views	Places the annotations for all views at the selected position.
Body Marker	---
Exam	Lists the available exams.
List	Displays the associated body markers.
All Bodymarkers	Displays all body markers.
Transducer Marker Location	Places the transducer marker at the selected position on the body marker.

Importing and Exporting Workflow Protocols

To import a protocol from a USB storage device:

Prerequisite: Connect an external storage device to the ultrasound system.

1. Access the configuration settings and then tap **Protocols**.
2. Click **Import**.
The system displays a list of protocols on the USB storage device and on the ultrasound system.
3. Select the protocols to import:
 - To select one or more protocols, select the check box for each required protocol.
 - To select all protocols, click **Select All**, if necessary.
 - To cancel the selection of all protocols, click **Unselect All**.
4. Click **Import Protocol**.

To export a protocol to a USB storage device:

Prerequisite: Connect an external storage device to the ultrasound system.

1. Access the configuration settings and then tap **Protocols**.
2. Click **Export**.
The system displays a list of protocols on the ultrasound system.
3. Select the protocols to export:
 - To select one or more protocols, select the check box for each required protocol.
 - To select all protocols, click **Select All**, if necessary.
 - To cancel the selection of all protocols, click **Unselect All**.
4. Select the required USB storage device and then click **Export Protocol**.

Contrast

Selection	Description
LVO Contrast	Provides the clip capture settings for left ventricle opacification contrast agent imaging.
Compression	Specifies the level of clip compression. ▪ Low compresses the frames of stored images with the low compression rate. Lower compression retains higher image quality. ▪ Medium compresses the frames of stored images with the medium compression rate. ▪ High compresses the frames of stored images with the high compression rate. Higher compression results in loss of image quality.
Capture Type	Specifies the clip store type. ▪ Prospective indicates a clip capture of the succeeding image. ▪ Retrospective indicates a clip capture of the previous image.
Clip Length	Provides measurement options for the duration of a clip capture. ▪ Time indicates the duration of each clip in seconds. ▪ Beat indicates the duration of each clip in heart beat cycles. ▪ R-Wave only indicates the capture point based on a delay beyond the R-wave.
GI Contrast	Provides the clip capture settings for general contrast agent imaging.
Compression	Specifies the level of clip compression. ▪ Low compresses the frames of stored images with the low compression rate. Lower compression retains higher image quality. ▪ Medium compresses the frames of stored images with the medium compression rate. ▪ High compresses the frames of stored images with the high compression rate. Higher compression results in loss of image quality.
Clip Length	Provides measurement options for the duration of a clip capture. ▪ Time indicates the duration of each clip in seconds. ▪ Beat indicates the duration of each clip in heart beat cycles.

Clip Store

Selection	Description
Compression	Specify the level of image compression and file size.
CaptureType	Specify the clip store type: ▪ Prospective indicates a clip capture of the succeeding image. ▪ Retrospective indicates a clip capture of the previous image.
Clip Length	Specify the length of the clip capture, measured in time or heart beats: ▪ Time indicates the time of each clip. ▪ Beat indicates the duration of each clip in heart beat cycles. ▪ R-Wave only indicates the capture point based on a delay beyond the R-wave.
Acquisition Rate	Defines the frame rate of the acquisition, for example, Normal or High .

Printer

Selection	Description
Printer	Note: To ensure images print correctly, select the name of the USB printer corresponding to each installed printer from the drop-down list.
PC Printer	Installed PC printers. Note: Printing a report or printing multiple images requires a universal PC printer or a black and white printer setting.
USB B&W Printer	Installed black and white printers.
USB Color Printer	Installed color printers.
Status	Status of the most recently assigned printer.
Thumbnails Printer	Output devices for printing thumbnails: ▪ PC Printer ▪ USB B&W Printer ▪ USB Color Printer ▪ None

Configuring and Installing Printers

Report printers are either directly connected to the ultrasound system or linked through a network path to the ultrasound system.

You must observe all safety precautions for connecting a non-medical device to the ultrasound system.

To add a network printer:

1. Access the configuration settings and then tap **Printer**.
2. Click **Install**.
The system displays the screen for installing the printer driver.
3. Click **Next**.
4. Click **Connect to a network printer** and then click **Next**.
 - To install a printer driver manually, click **Install a printer driver**, select a printer driver, and then click **Next**. The print drivers supported by the system only are displayed on the list.
 - The system automatically detects the available printers.
5. Click the required network printer and then click **Next**.
If the required printer is not detected, click **The printer that I want isn't listed**.
6. Follow the instructions on the screen.
7. Click **Close** when the installation is completed.
The system displays the message to restart the system.
8. Click **Yes**.
9. To confirm installation, select a printer and then print a report.
10. To save changes and exit the configuration menu, click **Save**.
11. To exit without saving your changes, click **Cancel**.

To select the PC printer or USB printer:

1. Access the configuration settings and then tap **Printer**.
2. Select the required printer from the drop-down list.

Configuring Printing Preferences

You can change the printing preferences, for example, paper type and size, page layout and orientation, and other settings available for your printer.

To change your print settings:

1. Access the configuration settings and then tap **Printer**.
2. Select the printer status settings.
A dialog box displays the printing preference selections.
3. Select the required options to optimize your print results.
4. Exit the dialog box and then save the configuration settings.

Barcode

Selection	Description
Barcode Settings	(Requires a connected barcode reader and your facility's representative barcode sample) Map the barcode sample to entry fields on the ultrasound system's patient registration form.
Scanned Data	Displays the characters transferred from the scanned barcode sample. The position of the cursor is listed beside the barcode characters. Note: You can position the cursor in the barcode to determine the starting position (offset) of each entry field.
Offset	Enter the starting position of the selected field in the barcode sample.
Length	Enter or modify the character length of the selected field. Note: The length field supports up to 65 characters.
Set	Confirm the offset and length values.
Sample Text	Enter or modify the sample text.
Use Separator	Enable or disable a separator for each barcode entry field. Enter a character in the text box. Note: You cannot use letters or numbers as separators.
Order	Specify the sequence of information. Note: Available only when Use Separator is activated.
(check box)	Enable or disable the selected entry field for use with the barcode reader. ▪ Last Name ▪ First Name ▪ Patient ID ▪ Performing MD ▪ Sonographer ▪ Accession No.
Offset	Enter or modify the starting position of the selected field in the barcode sample.
Length	Enter or modify the character length of the selected field. Note: The length field supports up to 65 characters.
Card Reader Setting	Specify the card reader settings to provide the direct entry of patient information in the registration form.
Scanned Data	Displays the characters transferred from the scanned card reader sample.
Patient ID	Displays the patient identification number transferred from the scanned card reader sample.
Offset	Enter the starting position of the selected field in the card reader sample.
Length	Enter the character length of the selected field. Note: The length field supports up to 65 characters.

Wired and Wireless Selections

Prerequisite: A working knowledge of networking principles is beneficial for configuring network connections.

The ultrasound system, or local host, can connect to a local area network (LAN) through a wired or a wireless connection. When both wired and wireless connections are configured and available, the ultrasound system gives priority to the wired connection.

Configuring Wired Connections to a Network

Consult your network administrator for the required configuration settings and ensure a network cable is connected to the network and to the ethernet port on the ultrasound system.

Configuring Wireless Connections to a Network

The ultrasound system uses network profiles (aliases) to store configurations for wireless network connections. You can enable and disable connections to wireless networks without a password even when changes to the wireless network options are protected by a password.

To configure a wireless connection to a network:

Prerequisite: Connect the wireless adapter to the ultrasound system before system startup.

1. Access the configuration settings and then tap **Wireless**.
2. To enable or disable connections to wireless networks, click **On** or **Off**.
3. To install a wireless adapter, click **Install...**
The system displays the device driver installation window.
 - a. Follow the instructions on the image screen.
 - b. Click **Close** when the installation is complete.
4. To add a network, enter the following network settings and then click **Apply**.

Note: Consult your network administrator for the required settings.

Selection	Description
IP Address Settings	---
DHCP	Enable the Dynamic Host Configuration Protocol (DHCP) if the ultrasound system uses dynamic IP addresses assigned by a DHCP server. The ultrasound system disables the IP Address , Subnet Mask , Default Gateway and MAC Address .
Static IP Address	Enable the Static Internet Protocol (IP) address if the ultrasound system uses the fixed IP address, and then enter the required settings for IP Address , Subnet Mask , Default Gateway .
Clear	Remove the entered IP address settings.
Apply	Apply the entered IP address settings.

To connect the ultrasound system to a network with a wireless connection:

1. Access the configuration settings and then tap **Wireless**.
2. Select the required wireless network from **Available Networks** and then click **Connect**.
3. To add a network profile, click **Add** and then enter the required settings.
4. To update the list of available networks, click **Refresh**.
The list displays the name of the network profile for the wireless network.
5. To modify a network profile:
 - a. Select the required wireless network from **Available Networks**.
 - b. Click **Edit** and then enter the required settings.
 - c. Click **OK**.

Network

You can configure the ultrasound system and connected devices, for example, servers for DICOM by creating, editing, and activating aliases (DICOM configurations).

See also: For information on configuring network settings, refer to Chapter 2 in this manual.

Utility

Selection	Description
Clip Playback	Enable or disable clip playback ▪ Play Only Active Clips disables simultaneous playback of all clips in the current study. ▪ Play All Clips enables simultaneous playback of all clips in the current study. ▪ Loop Aligned begins each playback cycles at the same time.
Slide Show	Configure the image duration and the number of replay loops ▪ Image Duration (sec) specifies the length of display for each image. ▪ Number of Replay Loops specifies the number of times each clip is played back.

Archive

Selection	Description
General	---
Options for Deleting Archived Studies	Options for managing the study history ▪ Delete Study History after Study is Deleted deletes the archived study history after the study is deleted. ▪ Manually Delete Study History after Study is Deleted provides the user options to delete the archived study history after the study is deleted. ▪ Keep Study History after Study is Deleted saves the archived study history after study is deleted.
Show Deleted Studies in Browser	Enable or disable the display of deleted studies in the patient browser
External Media	---
Copy Showcase to External Media when DICOM format checked.	Enable or disable automatic storage of the DICOM viewer software program with DICOM-formatted data archived to the external media.

Security

You can view the list of enabled and disabled security features on the ultrasound system. A check mark indicates an enabled feature, and an X indicates a disabled feature. The checklist also provides selections for configuring security settings.

Selection	Description
syngo Security	(Available only for the ultrasound system administrator) Indicates the status of the security feature ▪ Configure displays a dialog box to access user accounts. See also: For information about configuring user accounts, refer to Chapter 6 in this manual.
Windows 10	Indicates the operating system is Windows 10
Firewall	Indicates the firewall is configured
Ports	Indicates which encrypted ports are open or closed on the ultrasound system ▪ Detail displays a list of the names and descriptions of the ports. See also: Enabling Ports, page 1-26
Windows Hotfixes	Indicates which Windows hotfixes are installed on the ultrasound system ▪ Detail lists the name and date of the Windows hotfixes.
User Account Management	Indicates the ultrasound system includes user account management
Audit Trail	Indicates the ultrasound system is saving the security log files to the archive location. ▪ Detail displays additional selections for the audit trail. See also: For information about the audit trail for user accounts, refer to Chapter 6 in this manual.
FIPS (Federal Information Processing Standard)	(Available only for the ultrasound system administrator) Specifies the configuration of the Federal Information Processing Standard (FIPS) ▪ Enable enables or disables FIPS.
BIOS Password	Indicates the ultrasound system has a BIOS (basic input/output system) password
Disk Encryption	Indicates the data storage on the ultrasound system is encrypted
MDS2	Displays the Manufacturing Disclosure Statement for Medical Device Security ▪ Detail displays a file with details regarding the manufacturing disclosure statement for medical device security standard. – Print sends the file to the printer. – Export (USB) transfers a copy of the file to an external USB-compatible storage device.

Selection	Description
OSS Clearance	Displays the open source software included on the ultrasound system ▪ Detail displays a file with details regarding the open source software clearance. – Print sends the file to the printer. – Export (USB) transfers a copy of the file to the external USB-compatible storage device.
Client Certificate	(Available only for the ultrasound system administrator) Imports or exports a client certificate for wireless network connection. ▪ Configure displays a dialog box for exporting or importing the client certificate. See also: For information about importing and exporting the client certificate, refer to page 1-26.

Enabling Ports

(Available only for the ultrasound system administrator)

To open or close a port:

1. Access the configuration settings and then tap **Security**.
2. Click **Detail** for **Ports**.
3. Enable or disable a port.

Port Number	Port Group	Protocol
104	DICOM	TCP

Note: The following ports are enabled or disabled only by a Siemens Healthineers service representative.

Port Number	Port Group	Protocol
443	Remote Services	TCP

Importing or Exporting a Client Certificate

You can import or export a client certificate for use with wireless network.

To import or export a client certificate:

1. Access the configuration settings and then tap **Security**.
2. Click **Configure** for **Client Certificate**.
3. Follow the on-screen instructions provided in the Windows wizard to import or export the certificate.

Service

This selection provides the system serial number and software version number, along with password-protected service protocols and remote access for authorized Siemens Healthineers service representatives to troubleshoot and maintain the ultrasound system.

Selection	Description
Service	Displays the Siemens Healthineers Service Software window.
Remote	Displays the Remote Service Access Control window.
Restore Factory Settings	Deletes the user-defined settings and restores the factory default settings.

Note: The system displays text on the service screen in English only.

Enabling Anti-Virus Protection

You can enable or disable the anti-virus protection option.

See also: For information on enabling or disabling anti-virus protection, refer to Chapter 3 in the Instructions for Use.

Backing Up and Restoring Configuration Settings

You can backup and restore configuration settings and user-defined exam settings.

File Type Descriptions

Each file type saves a specific set of system settings.

Specify the ultrasound system files, for example, system configuration settings, protocols, or volume imaging parameters, to back up and restore.

To back up or restore configuration settings:

1. Access the configuration settings and then tap **Service**.
2. Click **Service** to access the settings in the service software and then click **OK** to display the user-configurable settings.
3. Click **Backup & Restore**.
4. Insert the external storage device into the USB port.

To	Do This
Save configuration settings and user-defined exam settings to a USB storage device	1. Click Backup from the Select list. 2. Select the check box for each required setting and then click Next. 3. Select the required storage device from the Destination list. 4. Select the check box for each required file type and then click Backup. The system displays a completion message when the process is finished. 5. Click OK. 6. Click X. Or, click Home.
Retrieve saved settings from USB storage device	Note: When you retrieve configuration settings, the system overwrites the existing configuration settings. When you retrieve a stored user-defined exam, the system overwrites any existing user-defined exam with the same name. 1. Click Restore from the Select list. 2. Select the required storage device from the Destination list. 3. Select the check box for each required file type and then click Restore. 4. Click OK to confirm the completion message. 5. To delete all existing user-defined exams on the system and add new user-defined exams: a. Select the Clean User-Defined Exam(s) on System check box in the Options section. b. Select the user-defined exams in the User-Defined Exam – Removable Media section. c. Click Restore>>. 6. To add user-defined exams to the system without deleting existing user-defined exams: a. Select the user-defined exams in the User-Defined Exam – Removable Media section. b. Click Restore>>. 7. Click Close. The system displays a message that power must be cycled to the system for completion of the operation. 8. Click OK to restart the ultrasound system.

2 Network Configuration

Configuring Network Settings	2
Prerequisites.....	2
Configuring a Local Host	3
Configuring a Network Export Server	5
Configuring a Destination Server.....	6
Verifying Network Settings	10

Configuring Network Settings

Use the configuration settings on the ultrasound system to configure the following devices:

Device	Description
Local Host	Local network settings for the ultrasound system.
Storage Server	Configuration settings for the devices where patient information, images and clips are stored. You can activate multiple storage servers and specify when images are saved to the server.
Worklist Server	Configuration settings for information system servers, for example, a hospital information system (HIS). You can configure the system to perform a search of procedures on the active worklist server.
MPPS Server	Configuration settings for the Modality Performed Procedure Step (MPPS) server that works with the worklist server to provide updates to exam status.
Printer Server	Configuration settings for multiple print servers.
Network Export Server	Configuration settings for the shared folder on the network.

The ultrasound system has integrated networked communication capability. Networked communication is in compliance with the Digital Imaging and Communications in Medicine (DICOM) 3.0 standard (also known as ACR/NEMA 3.0).

Prerequisites

- A working knowledge of Windows networking principles is necessary and beneficial for completing network configuration. An understanding of DICOM principles is necessary and beneficial for completing DICOM configuration.
- Close the patient study before performing configuration. Network configuration is disabled during a patient examination.
- Connect the ultrasound system, or local host, to a local area network (LAN) through a wired or a wireless connection.
- Network setup requires the following information for each configured device. Obtain these settings from your network administrator:
 - Ultrasound system host name, alias, Application Entity (AE) title, Internet Protocol (IP) address, subnet mask, and default gateway information is required
 - Media type and duplex mode required for connection to the host (ultrasound system)
 - Alias, application entity title, internet protocol address, and port number for each network device

Configuring a Local Host

You must configure local network settings for the ultrasound system (the local host) before configuring a network device or location.

Obtain the following information from your network administrator before configuring the local host:

- Alias
- Application Entity (AE) title
- Ultrasound system host name
- MAC address
- Internet Protocol (IP) address, not required if using DHCP
- Subnet mask, not required if using DHCP
- Default gateway information, not required if using DHCP

To configure local host settings:

1. Tap the configuration icon.



2. Tap **Network**.
3. Click **Local Host**.
4. To define a local host configuration, click **New**.
5. To edit an existing configuration, select the name (alias) from the list and then click **Edit**.
6. Configure the host type:
 - If the IP address is dynamic, or assigned by the server, click the **DHCP** check box.
 - If the IP address is static, click the **Static IP Address** check box.

7. Configure additional information, as required.

Selection	Description
IP Address	Internet Protocol (IP) address of the host (ultrasound system). Note: If DHCP is used for the export host, then reconfigure this field when the IP address changes.
Subnet Mask	Subnet mask for the host (ultrasound system). Typically identical to the subnet mask for the export host.
Default Gateway	Default gateway for the host (ultrasound system). Typically identical to the default gateway for the export host.
Alias	User-designated name for a host.
AE Title	Application Entity (AE) title for the configured device. Note: Disregard this entry. The system automatically populates this field if the ultrasound system is configured for DICOM. Otherwise, the system does not use this information.
Host Name	Name of local host.
MAC Address	Media Access Control (MAC) address detected for the host.
Show Network Status	Enable or disable the status of the network connection on the image screen.
Media Type	▪ Hardware Default ▪ Auto Select ▪ 10BASE-T ▪ 100BASE-TX ▪ 1000BASE-TX
Duplex Mode	▪ Hardware Default ▪ Half Duplex communicates in both directions, one direction at a time. ▪ Full Duplex communicates simultaneously in both directions.

8. To save changes, click **Apply**. To exit without saving the changes, click **Cancel**.
9. To remove an existing configuration, select the required name (alias) from the list and then click **Delete**.

Configuring a Network Export Server

You can create a shared folder on a destination computer, for example, a workstation, server, or personal computer.

To configure a network export location:

1. Tap the configuration icon.



2. Tap **Network**.
3. Click **Network Export Server**.
4. To define a network export location, click **New**.
5. To edit an existing configuration, select the name (alias) from the list and then click **Edit**.
6. Enter configuration settings for the location.

Note: The export host is the destination device.

Selection	Description
Auto Store to Network	Queues images for transfer to the shared folder on the network when you close the study.
Alias	System-defined "computer name" or user-defined name for the export host. This alias is not used by DICOM.
IP Address	Internet Protocol (IP) address of the configured device. Note: If DHCP is used for the server, then reconfigure this field when the IP address changes.
Ping Time Out in Seconds	Time limit in seconds for confirming a successful connection. Ping confirms a successful connection to the server. The system displays the results on the screen.
Shared Folder	Exact name of the shared folder on the export host, following standard Microsoft Windows conventions for specifying access.
Account	Account name used to log on to the export host. Examples: "archive1\jsmith" ("computer name" and user name for local account on export host) "jsmith" (user name for local account on export host) "domain1\jsmith" (domain and user name)
Password	Password used to log on to the export host, if applicable.
Test Write	Attempt storing "test" files to the export host.

7. To save changes, click **Apply**. To exit without saving the changes, click **Cancel**.
8. To remove an existing configuration, select the required name (alias) from the list and then click **Delete**.

Configuring a Destination Server

Obtain the following information from your network administrator before configuring the destination device (export host):

- Alias
- Application Entity (AE) title
- Internet Protocol (IP) address

To configure a destination server:

1. Tap the configuration icon.



2. Tap **Network**.
3. Click the required tab, for example **Storage Server**, **Worklist Server**, **Printer Server** or **MPPS Server**.
4. To define a server location, click **New**.
5. To edit an existing configuration, select the name (alias) from the list and then click **Edit**.
6. To remove an existing configuration, select the required alias from the list and then click **Delete**.
7. Enter configuration settings for the device.

Note: The export host is the destination device.

8. To configure storage server selections, configure additional information, as required:

Selection	Description
Store	Specifies the DICOM storage server.
SR	Specifies the structured report (SR) server to send a structured report.
Auto Store to DICOM	Queues images for transfer to the DICOM storage server: ▪ Store During Exam transfers each image when you store the image to the study. ▪ Store At End Of Exam transfers all images when you close the study.
Alias	System-defined "computer name" or user-defined name for the export host. This alias is not used by DICOM.
AE Title	Application Entity (AE) title for the configured device.
IP Address	Internet Protocol (IP) address of the configured device.
Port Number	Network port the configured device uses to communicate with the ultrasound system.
Time Out in Seconds	Time limit in seconds for confirming a successful connection. ▪ Ping confirms a successful connection to the server. The system displays the results on the screen. ▪ Echo confirms the application entity title and port number. The system displays the results on the screen.
Number of Times to Retry After Failure	Number of attempts to transfer data after an unsuccessful attempt.
Seconds Between Each Retry	Time interval in seconds after an unsuccessful attempt to transfer data before reattempting to transfer the data.
Write Time Out in Seconds	Time interval in seconds for sending data to the queue.
Connect Time Out in Seconds	Time interval in seconds for attempted connections between the host and the server.

Selection	Description
Image Format	Specifies the DICOM format used to transfer images to the storage server. Note: Calibration information is included only with the New Ultrasound format. ▪ Automatic specifies the highest quality format that is supported by the server: New Ultrasound, Old Ultrasound, or Secondary Capture (listed in sequence from highest quality to lowest quality). ▪ Original Image specifies the highest quality format that is supported by the server: Old Ultrasound (highest quality) or Secondary Capture (lowest quality). ▪ Secondary Capture specifies the Secondary Capture format. The Secondary Capture format does not support clips.
Monochrome	Specifies the monochrome conversion settings used to transfer the images to the storage server. ▪ No Monochrome conversion disables monochrome conversion for gray scale images and clips stored on the storage server. ▪ Store all Images and Clips as Monochrome stores all images and clips on the storage server in monochrome format. ▪ Store all Images as Monochrome stores all images on the storage server in monochrome format.
Storage Commit Server	Specifies when images are transferred to the storage server: ▪ After Each Image is Stored transfers each image to the storage server when you store the image to the study. ▪ After All Images are Stored transfers all images to the storage server when you close the study.

9. To configure worklist server selections, configure additional information, as required:

Selection	Description
Simple Search	(Applies to all configured worklist servers) When enabled, the system searches all procedures on the worklist server for the next 24 hours. Note: If you perform a streamlined search from the patient data form, the system uses data entered on the form to narrow the search.
Strict Verification	(Applies to all configured worklist servers) When enabled, the system verifies the DICOM standard tags on each study list retrieved from the worklist server.
Maximum number of worklist item	Maximum number of search results to display.

10. To configure MPPS server selections, configure additional information, as required:

Selection	Description
Store Image Format	Note: Calibration information is included with the New Ultrasound format only. ▪ New Ultrasound specifies the highest quality format that is supported by the server: New Ultrasound, Old Ultrasound, or Secondary Capture (listed in sequence from highest quality to lowest quality). ▪ Original image specifies the highest quality format that is supported by the server: Old Ultrasound (highest quality) or Secondary Capture (lowest quality). ▪ Secondary Capture specifies the Secondary Capture format. The Secondary Capture format does not support clips.

11. To configure the printer server selections, configure additional information, as required:

Selection	Description
Auto Print to DICOM	Queues images for transfer to the DICOM print server: ▪ Print when page is full prints an image when the page is full. ▪ Print at end exam prints all images when you close the study.
Orientation	Print orientation, for example, portrait or landscape
Display Format	Column and row format used for assembling pages to print
Copies	Number of copies to print
Film Size	File size to print
Printer Type	Printer type
Show advanced settings	Displays detailed settings for printer server
Medium Type	Specifies the print medium ▪ Paper selects paper. ▪ Clear film selects a clear film. ▪ Blue film selects a blue film.
Border Density	Defines the color of the region outside of the print area: typically black for film and white for paper
Min. Density	Note: Consult the print manufacturer for the appropriate value.
Max. Density	Note: Consult the print manufacturer for the appropriate value.
Print Priority	Print priority ▪ High indicates a high priority ▪ Medium indicates a medium priority ▪ Low indicates a low priority
Film Destination	Location of the film development ▪ Magazine indicates the exposed film sheet is developed in a darkroom. ▪ Processor indicates the exposed film sheet is fed directly into a processor and then developed.

Selection	Description
Empty Image Density	Density of the empty image area ▪ Black typically used for transparencies ▪ White typically used for paper
Polarity	Specifies the print color of the background and font. ▪ Normal prints black on white. ▪ Reverse prints white on black.
Magnification	Specifies how the system sends data to the camera if pixel interpolation is necessary. Pixel interpolation can be necessary when sending data to a high-resolution camera. ▪ Replicate indicates the interpolated value is identical to the last pixel value. ▪ Bilinear indicates the interpolated value is calculated linearly between the last pixel value and the next pixel value. ▪ Cubic indicates the interpolated value is calculated on a curve between the last pixel value and the next pixel value. This option is slowest but preserves fine detail better than bilinear interpolation. ▪ None indicates the image is left as-is.
Configuration	Note: This string is sent to the printer exactly as it displays in the field. This field is not required.

12. To save changes, click **Apply**. To exit without saving the changes, click **Cancel**.

Verifying Network Settings

You must test network configuration settings when you add a network server or edit an existing configuration.

To verify network settings for an alias:

1. Select the required alias.
2. To validate an IP address, click **Ping**.
Note: Your network administrator may have disabled the **Ping** function on your network. Even if **Ping** is disabled, the **Echo** test will still work.

The system displays the results of the test on the screen.

3. To validate the application entity and port number, click **Echo**.
Note: An Echo test is performed to verify that the device or server is properly configured for the ultrasound system.

The system displays the results of the test on the screen.

4. If the **Echo** verification fails:
 - a. Check the IP address, Application Entity Title, and port number.
 - b. Verify the application entity title of the ultrasound system is registered on the applicable server, if required.
5. Repeat steps 2 through 4 for each new or edited network configuration.

3 Measurement and Report Configuration

General Measurement Configuration	2
Measurement Display Configuration	5
Measurement Display for Obstetrics	6
Measurement Display for Cardiac	8
Measurement Display for Gynecology	8
Measurement Display for Pediatric	8
User-defined Measurement Selections	9
Touch Screen Selections	9
Measurement Labels	10
User-defined Measurement Labels	11
DICOM Structured Reporting	12
Obstetric Table	15
User-defined Formulas	16
User-defined Formulas for the Obstetric Exam	16
Reports	18
Report Section	18
Anatomy Assessments	18
Comments	19
Print Preview	19
Store Report	20

General Measurement Configuration

You can customize measurements for each exam type.

To configure the general measurement options:

1. Tap the configuration icon.



2. Tap **Measurement**.
3. Click **General**.
4. Enter the required settings.

Selection	Description
Guide Line	---
2D Caliper	Select the check box to display a dotted line between the anchored beginning and ending measurement markers.
Doppler / M-Mode	Select the display of lines indicating measurement guides. ▪ On displays a line. ▪ Off removes the line.
Marker Display Settings	Specify the measurement marker display settings.
Size	Select the measurement marker size, for example, small or large.
Type	Select the display of the measurement marker on the image screen: ▪ Object displays different shapes for each measurement marker ▪ Number displays consecutive numbers for each measurement marker
Caliper Default Position	Select the location on the image screen of the first measurement marker: ▪ Center displays the first marker in the center of the image screen but does not include a depth value. ▪ Depth displays the first marker in the center of the image screen, with a dotted line representing the depth from the skin line. The depth value displays on the image screen until you anchor the first marker.
Heart Rate Tool	Select the number of heart cycles.
Wall Motion Scoring	---
Segment Model	Select the segment model, for example, 16 segments or 17 segments.
Report View Type	Select the display of the report, for example, view plane, bull's eye, or both.
Measurement Results Display Settings	Select the default measurement method for each imaging mode.
2D, M, or Doppler (tab)	Select the required mode and specify the required display settings.
Default Position	Select the default location to display the measurement results, and position the pointer over the required location.
Background Opacity	Select the background transparency for measurement results display.

Selection	Description
Units and Precision (General Tools)	Specify the unit and precision settings for general measurements.
Default Tool	Specify the default measurement method by mode.
2D mode	Specify the default measurement method: ▪ Distance selects the default measurement method for measuring distance, for example, Open Trace ▪ Multi-point Distance selects the default measurement method for a measuring multiple distances, for example, 4 Point Dist ▪ Area selects the default measurement method for measuring an area, for example, Trace ▪ Circumference selects the default measurement method for measuring a circumference, for example, Ellipse ▪ Stenosis selects the default measurement method for measuring stenosis, for example, D-%Stenosis ▪ Angle selects the default measurement method for measuring an angle, for example, Cross Angle ▪ Volume selects the default measurement method for measuring a volume, for example, Vol 1 Dist to measure one volume distance ▪ Disk Volume selects the default measurement method for measuring a method of disk volume, for example, Disk Trace ▪ Ratio selects the default measurement method for determining a ratio, for example, Dist Ratio to measure a distance ratio ▪ Average selects the default measurement method for determining an average of multiple distance measurements, for example, 2 Distance Avg ▪ AP selects the default measurement method for measuring the pocket depth of amniotic fluid, for example, AP Depth
M-mode	Specify the default measurement method: ▪ Distance selects the default measurement method for measuring a distance, for example, 4 Point Dist ▪ HR selects the default measurement method for measuring a heart rate ▪ Slope selects the default measurement method for measuring a slope ▪ Time selects the default measurement method for measuring time ▪ Dist Ratio selects the default measurement method for determining a distance ratio ▪ Time Ratio selects the default measurement method for determining a time ratio

Selection	Description
Doppler	Specify the default measurement method: ▪ Velocity selects the default measurement method for measuring velocity, for example, PI Auto ▪ Acceleration selects the default measurement method for measuring acceleration or deceleration ▪ Velocity Ratio selects the default measurement method for determining a velocity ratio ▪ Time selects the default measurement method for measuring time ▪ Time Ratio selects the default measurement method for determining a time ratio ▪ RI selects the default measurement method for determining a resistive index ▪ A-Flow Volume selects the default measurement method for measuring a flow volume using an area ▪ D-Flow Volume selects the default measurement method for measuring a flow volume using a distance ▪ HR selects the default measurement method for measuring a heart rate ▪ VTI selects the default measurement method for measuring velocity time integral ▪ VTI Trace selects the default measurement method for measuring a velocity time integral trace. ▪ PHT Guide selects the default measurement method for measuring a valve area by pressure half time ▪ PHT selects the default measurement method for measuring a valve area on the irregular spectrum by pressure half time ▪ E/A Vel + DT selects the default measurement method for measuring a pressure half time and deceleration time ▪ dP/dt selects the default measurement method for measuring rate of change of pressure over time.
A-%Stenosis	Specify the default measurement method for measuring stenosis, for example, Ellipse
A-Flow Volume	Specify the default measurement method for measuring a flow volume of an area, for example, Trace

5. Click **Close**.

Measurement Display Configuration

You can configure the display of information on the measurement menu and on the patient report.

For other exam types you can configure the measurement display options only.

Note: To prevent the exam-specific settings from being overridden by the general settings, configure the general settings first and then customize the exam-specific settings.

To configure the measurement display options:

1. Tap the configuration icon.



2. Tap **Measurement**.
3. Click **General** and then click **Exam Specific**.
4. Select the exam from the **Application** drop-down list.
5. Enter the required settings.

Selection	Description
Measurement Results Display Settings	Select the default measurement method for each imaging mode.
Default Position	Select the default location to display the measurement results, and position the pointer over the required location.
Background Opacity	Select the background transparency for measurement results display.
Units and Precision (Measurement Labels)	Specify the unit and precision settings for exam-specific measurements. ▪ Apply to All Applications applies the specified measurement unit and precision settings to all exams.

Measurement Display for Obstetrics

Selection	Description
Add 1 week to EDD	Add an additional week to the estimated date of delivery. Note: You can change this setting during an exam. The ultrasound system updates the estimated date of delivery, gestational age, or last menstrual period in the report and patient data based on the changed setting. The system does not update information in structured reports previously sent to a network or the patient data stored on images or clips.
Growth Authors	Specify the settings for displaying the obstetric calculations in the measured results and growth analysis graphs for all reference authors. ▪ Percentile displays the percentile for the obstetric calculations. ▪ SD displays the standard deviation for the obstetric calculations.
OB	---
Growth Analysis Graphs	---
Link Studies	Select the segment model, for example, Manual or Automatic .
Graph Layout	Select the graph layout format. ▪ 1x1 displays one graph on a page in the patient report. ▪ 2x2 displays four graphs on a page in the patient report.
HC/AC Measurement	Select the measurement method for the head circumference and abdominal circumference ratio. ▪ HC with BPD/OFD measurement displays the values for biparietal diameter and occipital frontal diameter from the head circumference measurement. ▪ AC with TAD/APAD measurement displays the values for transverse abdominal diameter and abdominal anteroposterior diameter from the abdominal circumference measurement.

Selection	Description
Display SD / Growth Percentile	---
Notice of estimated Age	Specify the upper limit of the range for the standard deviation.
Display GA	Display or hide the gestational age.
Display EDD	Display or hide the estimated date of delivery.
Auto OB	Enable or disable automatic obstetric measurement.
BPD Outer to Outer	Select the check box to enable the biparietal diameter outer-to-outer measurement.
Early OB	---
Growth Analysis Graphs	---
Link Studies	Select the segment model, for example, Manual or Automatic .
Graph Layout	Select the graph layout format. ▪ 1x1 displays one graph on a page in the patient report. ▪ 2x2 displays four graphs on a page in the patient report.
Display SD / Growth Percentile	---
Notice of estimated Age	Specify the upper limit of the range for the standard deviation.
Display GA	Display or hide the gestational age.
Display EDD	Display or hide the estimated date of delivery.
Auto OB	Enable or disable automatic obstetric measurement.
BPD Outer to Outer	Select the check box to enable the biparietal diameter outer-to-outer measurement.
OB(Adv)	---
Growth Analysis Graphs	---
Link Studies	Select the segment model, for example, Manual or Automatic .
Graph Layout	Select the graph layout format. ▪ 1x1 displays one graph on a page in the patient report. ▪ 2x2 displays four graphs on a page in the patient report.
Display SD / Growth Percentile	---
Notice of estimated Age	Specify the upper limit of the range for the standard deviation.
Display GA	Display or hide the gestational age.
Display EDD	Display or hide the estimated date of delivery.
Auto OB	Enable or disable automatic obstetric measurement.
BPD Outer to Outer	Select the check box to enable the biparietal diameter outer-to-outer measurement.
Multi-Print	---
Configure Multiple Print Format	Specify the multiple print format settings.

Measurement Display for Cardiac

Selection	Description
Cardiac / Ped Echo / Neo Echo / TEE	---
General	---
HR Measurement Method	Select the measurement method used to detect heart rate: ▪ ECG ▪ Doppler / M-Mode ▪ Auto Trace
Default RAP	Select default value for right atrial pressure.
eSie Measure Setting	---
eSie Measure	Enable or disable the display of eSie Measure.
Display of Heart Cycles	Select the number of heart cycles for display in each imaging mode.

Measurement Display for Gynecology

Selection	Description
GYN	---
Auto Follicle Config	Specify the automatic follicle measurement settings.
Auto Follicle	Enable or disable the automatic follicle measurement.

Measurement Display for Pediatric

Selection	Description
Pediatric Hip GRAF Setting	---
Sonometer on screen option	Show or hide the GRAF sonometer.

User-defined Measurement Selections

You can configure the display of selections available during the measurement function.

To configure measurement function selections:

1. Access the configuration settings and then tap **Measurement**.
2. To create user-defined selections, click **Add New Style**.
3. To edit user-defined selections, click **Edit Style Name**.
4. To delete user-defined selections, click **Delete Style**.
5. Select the required exam type and user-defined selection, and then configure the following selections:

Selection	Description
Touch Screen	Assign selections to the controls corresponding to touch screen selections.
Left Panel	Define the display of settings for measurements, calculations, measurement tools, and measurement labels by exam type and operating mode.
Edit Label	Define the customized measurement and calculation labels by exam type and operating mode.
Formula	Define the customized measurement formula by exam type.
OB Table	(Available for all obstetric exam types) Configure the author reference settings for obstetric measurements.

Touch Screen Selections

You can configure the display of the controls corresponding to the touch screen selections.

To assign the touch screen selections:

1. Select the **Touch Screen** tab.
2. Select the operating mode.
3. Configure the following selections:

To	Do This
Assign the controls corresponding to the touch screen selections	1. Select a measurement type from the list and then select a label. 2. Click the right or left arrow.
Create a measurement group	1. Click Add Group. 2. Select a label from the list and then click the right or left arrow. 3. To edit a measurement group, click Edit Group. 4. To delete a measurement group, click Delete Group.
Change the location of two selections	▪ Click the selections and then click Swap.
Restore the default settings	▪ Click Restore Defaults.
Enable the system to advance to the next measurement label after you complete a measurement	1. Select the Measure Continuously check box. 2. Select the check box for each required measurement label in the required sequence. The system displays the selected measurement in the same color.

Measurement Labels

You can configure the location and sequence of measurement labels in the left panel of the image screen.

To configure the list of measurement labels:

1. Select the **Left Panel** tab.
2. Select the operating mode.

To	Do This
Assign the selections in the left panel	1. Select a measurement type from the list and then select a label. 2. Click the right or left arrow.
Create a measurement group	1. Click Add Group. 2. Select a label from the list and then click the right or left arrow. 3. To edit a measurement group, click Edit Group. 4. To delete a measurement group, click Delete Group.
Change the display order	▪ Select the label and then click the up or down arrow.
Display the measurement results on the image screen	▪ Select a label and then click the right arrow.
Restore the default settings	▪ Click Restore Defaults.
Enable the system to advance to the next measurement label after you complete a measurement	1. Select the Measure Continuously check box. 2. Select the check box for each required measurement label in the required sequence. The system displays the selected measurement in the same color.

User-defined Measurement Labels

You can create user-defined measurement labels.

To create a user-defined label:

1. Select the **Edit Label** tab.
2. Select the exam type.
3. Click **New**.
4. Specify the following information and then click **Save**.

Selection	Description
Label	Enter the name.
Exam	Select the exam.
Fetal	(Available only for obstetric, advanced obstetric, early obstetric, and fetal echo exams) When enabled, the ultrasound system generates the obstetrical measurement label for all fetuses.
Lt./Rt.	Select the right or left side of the anatomy for the measurement results.
Prox/Mid/Dist	Select the proximate, middle, or distant view for the measurement results.
Data shared among exams	(Available only with multiple exams) Select the option to display the measurement results of multiple exams.
Report Section	Specify the section and related measurement label to display in the report.
Highlight Value	Select the option to highlight the measurement results on the image screen.
Display Area	(Available only with calculation labels) Select the measurement result display area or report.
Value Display Type	Select the measurement result display type.
Supporting Tool	---
2D Mode / M mode / D mode (option)	Select the operating mode.
Available Method	Select the measurement method.
Name / Units / Precision (list)	Define the name, units, and precision for the selected method.

DICOM Structured Reporting

You must map a user-defined measurement label to include the measured result in DICOM structured reporting (SR). The ultrasound system provides sections and codes for mapping a measurement label to a DICOM SR template structure.

You can edit the SR mapping of system-defined measurement labels.

To include a measurement label in DICOM structured reporting:

1. Select the **Edit Label** tab.
2. Click the exam type.
3. Click the required operating mode tab: **2D**, **M**, or **Doppler**.
4. Select a user-defined measurement label from the list.
5. To map the measurement label to the SR structured report, click **SR Mapping**.
6. To edit the information for the measurement label and then map the label to the SR structured report:
 - a. Click **Edit**.
 - b. Modify the information for the measurement label and then click **Save and SR Mapping**.
7. Select the structured report template for the measurement label.

Selection	Description
SR Template	Specifies the DICOM structured reporting (SR) template.
(SR Template)	Identifies the DICOM structured reporting (SR) template. ▪ OB-GYN selects template identification (TID) 5000 OB-GYN Ultrasound Procedure Report. ▪ Fetal Echo selects TID 5220 Pediatric, Fetal and Congenital Cardiac Ultrasound Reports. ▪ Vascular selects TID 5100 Vascular Ultrasound Procedure Report. ▪ Cardiac selects TID 5200 Cardiac Ultrasound Procedure Report.
(SR Template Section)	Lists available sections for the selected SR template. ▪ OB-GYN SR Template Section ▪ Fetal Echo SR Template Section ▪ Vascular SR Template Section ▪ Cardiac SR Template Section Note: If the template identification (TID) sections are not applicable to the measurement label, a user-defined results section is available.
Anatomy or Vessel	(Available only for double concept codes) Lists predefined anatomical structure or a vessel for the structured report. The predefined anatomy or vessels are based on the selected SR template section.

Selection	Description
Coding Scheme Designator	Identifies the source of the code. ▪ LN indicates the mapping is defined by Logical Observation Identifiers Names and Code (LOINC). ▪ SRT indicates the mapping is defined by Systematized Nomenclature of Medicine Reference Terminology (SNOMED RT). ▪ DCM indicates the mapping is defined by Digital Imaging and Communications in Medicine (DICOM). ▪ 99SIEMENS indicates the mapping is defined by the factory. ▪ Use 99 as a prefix to indicate a private coding scheme designator, for example, 99NAME, where NAME is the name of your facility.
Add	(Available only when configuring a user-defined SR template) Adds a user-defined coding scheme designator.
Delete	(Available only when configuring a user-defined SR template) Removes the user-defined coding scheme designator.
Code Value	System-generated based on the selected predefined anatomy or vessel. Note: If you created a user-defined anatomy or vessel, you must also create a user-defined code value. The maximum number of characters is 16.
SR Tool Setting	Specifies the tool result for the measurement label.
Tools & Outputs	Lists the measurement tools assigned to the measurement label and each tool result. Note: Map each tool result for structured reporting.
Output	Displays the selected tool result.
Text String for SR	Adds a user-defined description of the tool result.
Concept Code Meaning	Provides a list of concept code meanings for a tool result. When you select a predefined meaning, the system completes the following fields: Coding Scheme Designator, Code Value. Note: If you create a user-defined concept code meaning, you must also specify the coding scheme designator and code value.
Add	(Available only when configuring a user-defined concept code meaning) Adds a user-defined coding scheme designator.
Delete	(Available only when configuring a user-defined concept code meaning) Removes the user-defined coding scheme designator.

Selection	Description
Code Value	System-generated based on the selected predefined anatomy or vessel. Note: If you created a user-defined anatomy or vessel, you must also create a user-defined code value. The maximum number of characters is 16.
Measurement Orientation	Specifies the position of the measurement, for example, longitudinal.
Postprandial	(Available for the vascular SR template) Specifies the measurement is taken after a meal.
Measurement Methods	(Available for the fetal echo and cardiac SR templates) Specifies the echo measurement method, for example, Planimetry.
Finding Site	(Available for the fetal echo and cardiac SR templates) Specifies the echo anatomic site, for example, Cardiac valve annulus.
Flow Direction	(Available for the fetal echo and cardiac SR templates) Specifies the direction of blood flow.
Cardiac Cycle Point	(Available for the fetal echo and cardiac SR templates) Specifies the cardiac phase.
Image View	(Available for the fetal echo and cardiac SR templates) Specifies the echocardiography image view.

8. To exit while SR mapping is in progress, click **Cancel**.
The ultrasound system saves only the user-defined measurement label.
9. Click **Finish**.
The ultrasound system saves the measurement label and the SR mapping definitions.

Obstetric Table

(Available for all obstetric exam types)

You can customize the list of reference authors for obstetric measurements and calculations.

To select a reference author:

1. Select **OB**.
2. Select the **OB Table** tab.
3. Specify the following information:

Selection	Description
Gestational Age	Select the author references for gestational age.
Growth Analysis Graphs	Select the author references for the growth analysis graph.
Ratio	Select the author references for ratio.
Estimated Fetal Weight	Select the author references for estimated fetus growth.
Same Author for All Labels	Select the author reference and apply the same author reference to gestational age.

4. Click **Close**.

The system applies the reference authors for obstetric measurements and calculations to all obstetric exams.

User-defined Formulas

You can create user-defined formulas.

To create a user-defined formula:

1. Access the configuration settings and then tap **Measurement**.
2. Click the exam type.
3. Click **New** and then enter the required information.
4. To edit a user-added formula, select the required formula and then click **Edit**.
5. To delete a user-added formula, select the required formula and then click **Delete**.
6. To make a correction in the formula, position the cursor after the incorrect element and then press the backspace key on the keyboard.

Selection	Description
Custom (tab)	---
Reference	Specify the reference settings for obstetric measurement result.
Formula Name	Enter the label name.
Output Value	Specify the output value including unit and precision.
Formula	Define the measurement formula: ▪ Constant specifies the labels and corresponding value. ▪ Variable specifies measurement labels, units and prefix by the imaging mode and exam type. ▪ Operation specifies arithmetic operation for calculation.

User-defined Formulas for the Obstetric Exam

The system displays the formula in the text box next to the formula name as you enter each element. Example: 'A*LOG('BPD'+COS('ATD'))-1.8327.

You can use parentheses, but do not use spaces to separate elements in the formula. Do not delete quotation marks entered by the system. You can enter up to 64 characters. To conserve space, you can define a numerical value as a constant.

To create a user-defined formula for obstetrics:

1. Select **OB**.
2. Click the **Formula** tab and select from the following formula tabs:
 - **Gestational Age**
 - **EFW**
 - **Growth Analysis Graphs**
 - **Ratio**
 - **Custom**
3. To create a formula, click **New**.
4. To edit a user-defined formula:
 - a. Select the required formula and then click **Edit**.
 - b. To make a correction in the formula, position the cursor after the incorrect element and then press the backspace key on the keyboard.
5. Specify the following information:

Selection	Description
Reference	Specify the reference settings for obstetric measurement result for the selected formula: ▪ Reference (Author) selects the reference author for determining the measurement results. ▪ Variable Label selects the variable label. ▪ Range specifies the lower and upper limits in weeks and days of the selected variable used to calculate the fetal growth. ▪ Data Type selects the layout for the data, for example, Table or Formula. ▪ Result Unit selects the units for the measurement results, for example, Week or Day. ▪ Displayed Value selects the unit and precision. ▪ Label Order selects the sequence for the measurement labels. ▪ Precision specifies the precision for the measurement value. ▪ Deviation selects the range for the fetal growth for the fetal growth pattern. ▪ Label Name inserts the name of the label. ▪ Output Value specifies the output value for the unit and precision. ▪ Measurement Result selects the option to calculate the fetal growth for the selected measurement label and reference author.
Formula	Define the measurement formula: ▪ Constant specifies the labels and corresponding value. ▪ Variable specifies measurement labels, units and prefix by the imaging mode and exam type. ▪ Operation specifies arithmetic operation for calculation.
SD Limits Formula	Define the standard deviation limit formula: ▪ Constant specifies the labels and corresponding value. ▪ Variable specifies measurement labels, units and prefix by the imaging mode and exam type. ▪ Operation specifies arithmetic operation for calculation.

6. To delete a user-defined formula, select the formula and then click the delete button.
7. Click **Save**.

Reports

You can configure patient report settings and insert user-defined information.

Report Section

You can show or hide sections in the patient report.

To configure a report section:

1. Access the configuration settings and then tap **Measurement**.
2. Click **Report**.
3. Click **Report Section**.
4. Click the exam type.
5. To display a section in the report, select the report section check box.
6. To reposition a section in the report, select the report section and then click the up or down arrow.
7. Click **Close**.

Anatomy Assessments

Anatomy assessments provide descriptions on the patient report, for example, for assessing the condition of an anatomical structure. An anatomy assessment is an area of interest, for example, a fetal position. A selection is a description of the area of interest, for example, breech.

The system identifies user-defined assessments with an asterisk (*).

To configure anatomy assessments:

1. Access the configuration settings and then tap **Measurement**.
2. Click **Report**.
3. Click **Anatomy Assessment**.
4. Click the exam type.

To	Do This
Create a user-defined anatomical assessment or selection	1. Enter a name. 2. To create a new description: a. Click Assessment Order. b. Enter a description. 3. To assign a description to an item: a. Click Selection Order. b. Enter a description. 4. Click Add.
Delete a user-defined anatomical assessment or selection	1. Select the area of interest or the description. 2. Click Delete. The item or selection is only removed from report.
Add user-defined anatomical assessment and selection	▪ Select the anatomical assessment check box.
Reposition an anatomical assessment in the list	1. Select the area of interest or the description. 2. Click the up arrow or down arrow.

Comments

You can create up to 20 comments for the patient report.

To configure a comment for the patient report:

1. Access the configuration settings and then tap **Measurement**.
2. Click **Report**.
3. Click **Comments**.
4. Select the exam type and then enter a comment.

Print Preview

You can define a logo and a signature, and customize the display of information for the institution and patient in the patient report.

To configure the print layout for the patient reports:

1. Access the configuration settings and then tap **Measurement**.
2. Click **Report**.
3. Click **Print Preview**.
4. Click the exam type and then specify the following information:

Selection	Description
Site Information	Specifies the logo for the selected exam. ▪ Load inserts an image. ▪ Remove deletes an image. ▪ Text 1 inserts text for display with the logo. ▪ Font specifies the font for the image.
Signature	Specifies a signature for the selected exam. ▪ Add inserts a signature. ▪ Remove deletes a signature.
Section Display	When enabled, hides information in the patient report. ▪ Hide institution hides institution name and physician- or sonographer-identifying information. ▪ Hide patient information hides patient name, identification number, age, gender, height, weight, body surface area (BSA), blood pressure (BP), and study date. ▪ Hide exam specific information hides exam-specific information, for example, last menstrual period or blood pressure. ▪ Hide all measurement results hides measured results. ▪ Hide assessment hides anatomy assessments. ▪ Hide comments hides the comment section.

Store Report

You remove identifying information when you export a patient report.

To remove data from an exported patient report:

1. Access the configuration settings and then tap **Measurement**.
2. Click **Report**.
3. Click **Store Report**.
4. Select the required options.
 - **Hide institution** hides data identifying the healthcare institution.
 - **Hide patient information** hides patient data.

4 Clinical Applications Configuration

Configuring Stress Echo Settings.....	2
Creating and Modifying a Protocol	2
Creating and Modifying a Protocol Stage	3
Creating and Modifying a Protocol View.....	4
Exporting and Importing Protocols	5

Configuring Stress Echo Settings

You can create, edit, and delete stress echo protocols. The changes take effect the next time you activate the stress echo feature.

Creating and Modifying a Protocol

When you create or edit a protocol, you can configure the stages and the number of clips acquired for views in the stage.

To create or modify a protocol:

1. Tap the configuration icon.



2. Tap **Stress Echo**.

To	Do This
Copy a protocol	1. Select a protocol and then click Copy. 2. Enter a name and then click OK.
Add a protocol	1. Select a protocol and then click Add. 2. Enter a name and then click OK.
Rename a protocol	1. Select a protocol and then click Rename. 2. Enter a name and then click OK.
Modify a protocol:	1. Select a protocol and specify the following fields in Protocol Settings. – Check ROI On to display the region of interesting area in stress echo setup screen. – Check Allow Selection of Multiple Clips per View to select multiple clips. 2. Select the required Review Settings option. – Select the required display layout in the stress echo review screen from the Display Format drop-down list. – Select the required section of the clip to play when you activate playback from the Default Playback Mode drop-down list, for example, Full R-R, Adjusted, Diastole Only or Systole Only. – Select the required playback speed in stress echo review screen from the Default Playback Speed drop-down list. – Select the required Default Clip Sorting option, for example By View or By Stage.
Delete a protocol	▪ Select a protocol and then click Delete.

Creating and Modifying a Protocol Stage

To	Do This
Copy a stage	1. Select a stage and then click Copy. 2. Enter a name and then select a unit from the drop-down list. 3. Click OK.
Add a stage	1. Select a stage and then click Add. 2. Enter a name and then select a unit from the drop-down list. 3. Click OK.
Rename a stage	1. Select a stage and then click Rename. 2. Enter a name and then select a unit from the drop-down list. 3. Click OK.
Modify a stage:	1. Specify the following fields in Clip Settings. – Select the required option from the Clip Store Type drop-down list: Retrospective, Prospective, or Continuous. – For non-continuous stage, select the number of clips stored when acquisition is activated from the Clips per Capture. The maximum number of clips for retrospective is 4, and for prospective is 5. – Select the type of acquisition in Beat or Time in Clip Length. – Select the type of capture offset in Time and R-Wave in Capture Offset. 2. Specify the following fields in Additional Settings: – To specify the time to enable Start Stage Timer, select At Start of Stage or At First Clip Store. – To disable Start Stage Timer, select No. – To specify the location to enter the review screen, select After View Capture or After Stage Capture in Enter Review Screen. – To disable Enter Review Screen, select None.
Delete a stage	▪ Select a stage and then click Delete.

Creating and Modifying a Protocol View

To	Do This
Copy a view	1. Select a view and then click Copy. 2. Select the required option. – Select Pre-defined List and then select the pre-defined view from the drop-down box. – Select Custom and then enter a name for the copied view. 3. Click OK.
Add a view	1. Select the view and then click Add. 2. Select the required option. – Select Pre-defined List and then select the pre-defined view from the drop-down box. – Select Custom and then enter a name for the copied view. 3. Click OK.
Rename a view	1. Select a view and then click Rename. 2. Select the required option. – Select Pre-defined List and then select the pre-defined view from the drop-down box. – Select Custom and then enter a name for the copied view. 3. Click OK.
Delete a view	▪ Select a view and then click Delete.

Exporting and Importing Protocols

To export protocols:

1. Connect an external storage device to the ultrasound system.
2. Click **Export**.
3. Select the detected external storage device and the required protocols and then click **Export**.

To import protocols:

1. Connect an external storage device to the system.
2. Click **Import**.
3. Select the detected external storage device and the required protocols and then click **Import**.

5 Export and Import

Exporting Patient Studies	2
Storing Images and Structured Report to Network Servers	3
Importing Patient Studies	4
Creating a Slide Show	4
Creating and Editing Teaching Files.....	5

Exporting Patient Studies

Use the patient browser to export closed patient studies to external media, for example, a USB storage device.

To export a patient study:

1. Press **PATIENT** and then tap **Browser**.
2. Select the time when the patient study was completed in the **Show Files** list and then select the location of the studies in the **Source Media** list, if necessary.
3. Specify the search options and then tap **Search**.
The system displays the patient information for available studies in the study list.
4. Select one or more studies in the study list and then tap **Export(Import)**.
5. Select the target location for the exported studies in the **Destination** list.
6. To remove patient-identifying information from the study, select the **De-identify** check box.
Note: You can enter text to replace the patient name in the text box.
7. Select the format for the exported images:
 - **DICOM** exports images in the DICOM (Digital Imaging and Communications in Medicine) format.
 - **Tiff/Avi** exports images in Tagged Image File Format and exports clips in Audio Video Interleaved format.
8. Select the options for disk media:
 - **Finalize** prevents additional storage to the inserted disk media.
 - **Eject** ejects the disk media.
9. Click **Export**.
The system displays a progress bar indicating the status of the export process.
10. When the export is complete, click **OK**.

Storing Images and Structured Report to Network Servers

You can send selected studies and included structured report to a connected server. You can also send stored images during the current study when the system is configured for in-progress store.

To send studies to a connected server:

Note: The selected studies must be stored on the hard disk of the ultrasound system.

1. Press **PATIENT**, and then tap **Browser**.
2. Select the time when the patient study was completed in the **Show Files** list, and then select the location of the studies in the **Source Media** list, if necessary.
3. Specify the search options, and then click **Search**.
The system displays the patient information for available studies in the study list.
4. Select one or more studies in the study list.
5. Select the target location in the **Output Device** list and then tap **Send**.
The system adds all images from the selected studies to the queue.
6. To verify the status of the queue, tap **Queue**.

See also: For information about the indication of the queue status, refer to Appendix B in the Instructions for Use manual.

Importing Patient Studies

Use the patient browser to import patient studies from external media, for example, a USB storage device.

To import a patient study:

1. Press **PATIENT**, and then tap **Browser**.
2. Select the location of the studies in the **Source Media** list.
3. Select one or more studies from the external media.
4. Click **Export(Import)**.
The system copies the study to the system hard drive and adds the patient information to the study list.
5. To confirm the studies were imported, select the hard drive in the **Source Media** list.

Creating a Slide Show

You can temporarily hide images on the review screen and display only selected images. You can also view selected images sequentially as a slide show.

Use the configuration settings to specify the length of time to display each image and the number of times to play a clip.

To create a slide show:

1. Choose a method to activate review:
 - To activate review during an exam, press **REVIEW**.
 - To activate review of a saved exam, activate the patient browser and then double-click the patient name in the study list.
2. Select the required images in the review screen.
3. To view selected images and clips as a slide show, tap **Slide Show**.
The system displays each image or clip in a full-screen display format.
4. To pause the slide show, click **Stop**.
5. To restart the slide show, click **Start**.
6. To stop the slide show and redisplay the previous screen, click **Exit**.

Creating and Editing Teaching Files

You can create a teaching file with images from multiple patient studies in a single file. The system removes all patient-identifying information from the images.

To create a teaching file:

1. Choose a method to activate review:
 - To activate review during an exam, press **REVIEW**.
 - To activate review of a saved exam, activate the patient browser and then double-click the patient name in the study list.
2. Select the required images in the review screen or thumbnail panel.
3. Tap **Teaching File**.
4. Enter a name for the teaching file.
5. Click **New**.

The system creates a new study with the selected images and adds the study to the patient browser.

To add images to a teaching file:

1. Choose a method to activate review:
 - To activate review during an exam, press **REVIEW**.
 - To activate review of a saved exam, activate the patient browser and then double-click the patient name in the study list.
2. Select the required images in the review screen or thumbnail panel.
3. Tap **Teaching File**.
4. Select the required file from the list and then click **Append**.
5. To add images from another study to the teaching file, open the required study and then repeat steps 2 through 4.

6 Security Configuration

Overview.....	2
Terminology Used in this Chapter	2
Configuring Security Settings	3
Categories of Settings for Security Configuration.....	4
User Account Management	6
Groups and Roles.....	6
Assigning Functional Privileges to Users.....	7
Adding a User Account.....	9
Editing a User Account.....	10
Configuring the User Account Password.....	10
Viewing Security Settings.....	11
Audit Trail Management	12
Configuring the Audit Trail	12
Archiving the Audit Trail.....	13

Overview

The information in this chapter is intended for the ultrasound system administrator.

Terminology Used in this Chapter

Term	Definition
User account	A login account assigned to each user to log in to the ultrasound system. User accounts are configured to restrict access to functional privileges on the ultrasound system.
meduser	A default internal user account.
Privilege	A restriction of access to functions or tasks on the ultrasound system; for example, sending data. User accounts can be assigned to roles to manage privileges. Privileges are also referred to as Functional Privileges .
User	A single user, identified on the system with a user account.
Group	Note: The system provides the group classification only. Permissions cannot be assigned to groups. A classification for a user account. Groups include user groups. The following list includes the default groups: ▪ Administrator/SecurityAdmins ▪ EmergencyAccess ▪ Standard
Role	A classification for a user account, providing access to functions or tasks on the ultrasound system. The following list includes the default groups and roles, which cannot be deleted: ▪ Administrator/SecurityAdmins provides accessibility based on the functional privileges and access to configure the ultrasound system. This group and role includes the ultrasound system administrator. ▪ EmergencyAccess provides access without a password for emergency login. The accounts have limited access. ▪ Standard provides access to all patient groups except service. The standard group and role cannot be deleted.
Audit trail	A collection of security log files. The security log file tracks user actions, including powering on the ultrasound system, logging in, accessing files, modifying patient data and images, and powering off the ultrasound system. Only authorized users can view the security log files. Security log files are saved to the system hard disk, a storage device, or a network location. Saving the security log file is also referred to as archiving the audit trail.
Federal Information Processing Standard (FIPS)	The Federal Information Processing Standard is a United States government computer security standard used to accredit cryptographic modules.

Configuring Security Settings

You can manage access to patient data on the ultrasound system by configuring user accounts.

Use the configuration settings on the ultrasound system to set up user accounts and passwords, configure groups and roles, update passwords, and configure security log files.

Categories of Settings for Security Configuration

The security package provides configuration options for managing user accounts and providing functional privileges.

User Management includes options for configuring users, groups, and roles.

- **Users** lists the individual user accounts. Use this selection to add and configure additional user accounts.
- **DICOM Nodes** lists the virtual accounts for transferring network data.
- **Internal Users** lists the internal accounts with access to all patient data and access to configure the ultrasound system.
- **Groups** lists the current group accounts. Use this selection to add and configure additional groups.
- **Roles** lists the current roles. Use this selection to add and configure additional roles.

Security Management provides options for assigning access to patient data and tasks on the ultrasound system.

- **Functional Privileges** lists access to functions and tasks on the ultrasound system.

To configure security settings:

1. Tap the configuration icon.



2. Tap **Security**.
3. Click the **Configure** selection for **syngo Security**.
4. Click the menu option for the required settings.
5. Enter the required settings.
6. Save the changes:
 - To save the settings, click **Apply**. Or, click **OK**.
 - To discard the changes and exit the configuration settings, click **Cancel**.

To enable the Federal Information Processing Standard:

1. Access the security configuration settings.
2. Click the **Enable** selection for the **FIPS (Federal Information Processing Standard)**.

User Account Management

You can manage the following user account configuration:

- Creating, editing, and deleting user accounts, including user names and passwords
- Configuring access to patient data and ultrasound tasks
- Integrating a group from a network domain
- Configuring, viewing, and archiving the security log file

Groups and Roles

You can create groups and roles to use in addition to the default groups and roles. Each name of the group or role must be unique and cannot include special characters. If a user group name and role name are identical, the system displays a message indicating all changes are applicable for both the user group and the role. After you create a group or a role, you cannot change the name.

To create or delete a user group or role:

1. Access the security configuration settings.
2. Click the **Configure** selection for **syngo Security**.
3. Click **Security Configuration**.
4. Select **User Management**.
5. To create a user group:
 - a. Click **Groups** and then click **More Actions**.
 - b. Select **New** and then select **Group**.
 - c. Enter a name and description.
 - d. Select the **Members** tab.
 - e. To add users, select the users from the **Available Users** list and then click **Add**.
 - f. To delete users, select the users from the **Assigned Users** list and then click **Remove**.
6. To create a role:
 - a. Click **Roles** and then click **More Actions**.
 - b. Select **New** and then select **Role**.
 - c. Enter a name and description.
 - d. Select the **Owners** tab.
 - e. Select a domain from the **Domain** drop-down list.
 - f. To add users, select the users from the **Available Users** list and then click **Add**.
 - g. To delete users, select the users from the **Assigned Users** list and then click **Remove**.
7. To delete a user-defined group or role:
 - a. Click the group or role and then click **More Actions**.
 - b. Select **Delete** and then select **Yes**.
8. Save the changes.

Assigning Functional Privileges to Users

You can assign functional privileges to a user account. The system displays an asterisk (*) next to any field with unsaved changes.

You can limit the list of groups and roles by selecting a filter at the top of the screen.

- A classification of user accounts includes the following selections:
 - **Everyone** includes all users and all groups
 - **Groups** includes a list of groups
 - **Roles** includes a list of roles
 - **Users** includes a list of individual user accounts
- The system displays a drop-down list of groups or roles based on the filter.

To assign functional privileges to users:

1. Access the security configuration settings.
2. Click the **Configure** selection for **syngo Security**.
3. Click **Security Configuration**.
4. Assign access to tasks on the ultrasound system by selecting the required options.
 - a. Select **Security Management** and then select **Functional Privileges**.
 - b. Select **By User / Role**.
 - c. Select the **Grant** check box next to the required selection and then enter or modify the configuration settings.

The access is specific to each user or role.

Selection	Description
Privilege	Lists the functions or tasks available to assign to the user account.
Effective	When the check box is selected, Indicates the effective access for the user account.
ArchNet (tab)	Export indicates access to export patient data. Send indicates access to send patient data.
AuditTrail (tab)	Archive indicates access to save the security log file. ConfigureReadWrite indicates access to configure, read, and edit the security log file. SetFilter indicates access to configure the security log file. View indicates access to view the security log file.
Common (tab)	RemoteAssistance indicates access to launch a remote service session. ServiceAccess indicates access to the local service page, including access to the backup and restore feature.

Selection	Description
CorrectRearrange (tab)	Correct indicates access to edit the patient registration information.
Filming (tab)	Expose indicates access to print patient data.
PatientBrowser (tab)	Delete indicates access to delete patient data.
PatientRegistration (tab)	InvokeRegistration indicates access to register a patient.
SecuritySystem (tab)	BypassBlockedSystem indicates access to bypass a system blocked as a result of the security log file filling the system hard disk. ConfigurationReadWrite indicates access to view and edit the security configuration. InteractiveLogin indicates access to log in when the system is locked. OverrideLoginErrors indicates a user can access the system when the system displays login errors; for example when an exam is in progress. SetMode indicates a user can access the system.

5. Save the changes.

Adding a User Account

You can add and configure user accounts, and change a user account password.

To add a user account:

1. Access the security configuration settings.
2. Click the **Configure** selection for **syngo Security**.
3. Click **Security Configuration**.
4. Create a user account.
 - a. Select **User Management**.
 - b. Click **Users** and then click **More Actions**.
 - c. Select **New** and then select **User**.
 - d. Enter the information in the **User** tab.

Selection	Description
Name	Login name for the user account. The name is not editable after you save the user account.
Full Name	Full name including first name and last name of the user
Description	Description of the user
New password	Password for the user
Confirm password	Password for the user
Account is disabled	When the check box is selected, the user account is not active and the user cannot log in using a user name and password.
Password never expires	When the check box is selected, the password for the user account does not expire.
User cannot change password	When the check box is selected, the user cannot modify the password for the user account.
Account is locked out	When the check box is selected, the user account is locked and the user cannot log in to the ultrasound system.

5. To add or remove the user to a group:
 - a. Select the **Member of** tab.
 - b. Select a group from the **Available Groups** list and then click **Add**.
 - c. To remove the user from a group, select the group from the **Assigned Groups** list and then click **Remove**.
6. To add or remove the user to a role:
 - a. Select the **Owner Of** tab.
 - b. Select a group from the **Available Roles** list and then click **Add**.
 - c. To remove the user from a role, select the role from the **Assigned Roles** list and then click **Remove**.
7. To delete a user account:
 - a. Click the user account and then click **More Actions**.
 - b. Select **Delete** and then select **Yes**.
8. Save the changes.

Editing a User Account

To edit a user account:

1. Access the security configuration settings.
2. Click the **Configure** selection for **syngo Security**.
3. Click **Security Configuration**.
4. Select **User Management**, select **Users** and then select the user from the list.
5. To filter the list of users:
 - Select a classification from the drop-down menu.
 - Enter the name of the new filter in the text box and then press enter on the keyboard.
6. To change a user password, select a user and then select the **User** tab.
7. Update the required information.
8. Save the changes.

Configuring the User Account Password

You can configure the password requirement for all user accounts, and reset a specific user account password.

To configure user account passwords:

1. Access the security configuration settings.
2. Click the **Configure** selection for **syngo Security**.
3. Click **Password Configuration**.
4. Select the required tab and then enter or modify the configuration settings.

Selection	Description
Reset Password (tab)	Displays a list of user accounts to reset the password.
Select All	Selects all of the user accounts on the list.
Clear All	Cancels the selection of all user accounts on the list.
Reset Password	Resets the password for the selected user accounts. The system displays a dialog box requesting confirmation to reset the password for the selected user accounts and restart the ultrasound system.
Password Complexity (tab)	Displays the password configuration requirements; for example, a minimum number of uppercase characters.
Default Settings	Selects the default setting for the password of the selected user accounts.
Account Policy (tab)	Displays account configuration requirements; for example, the duration to lock an account.

5. To save the settings, click **OK**.
6. To restart the ultrasound system, click **OK**.

Viewing Security Settings

You can view groups and roles with assigned functional privileges.

To view the security settings:

1. Access the security configuration settings.
2. Click the **Configure** selection for **syngo Security**.
3. Click **Security View**.
 - **Groups/Roles** shows the user groups and roles.
 - **Privileges** shows the functional privileges assigned to a user.

Audit Trail Management

Configuring the Audit Trail

You can configure the audit trail by selecting the events recorded in the security log file and specifying the archive location.

To configure the audit trail:

1. Access the security configuration settings.
2. Click the **Detail** selection for **Audit Trail**.
3. Click **Audit Trail Settings**.
4. Select the required destination for the security log file and then enter or modify the configuration settings.
5. To view the next configuration screen, click **Next**.
6. To view the previous configuration screen, click **Back**.

Selection	Description
Local file system	Selects the ultrasound system as the destination for the security log file.
Audit Trail directory	Selects the location for the archive of the security log file.
Max Size (kBytes)	Specifies the maximum size of the security log file.
Zipped	When enabled, compresses the security log file.
Disk quota checks	Warning level [%] specifies the storage capacity level of the security log file as a percentage of the maximum size, when enabled. Quota limit [%] specifies the quota limit of the security log file as a percentage of the maximum size, when enabled.
Total diskspace of Audit Trail drive	Displays the quota limit for audit trail in megabytes.
Used diskspace of Audit Trail drive	Displays the amount of the used disk space for audit trail in megabytes.
Central syslog server	Selects a network location as the destination for the security log file.
Host	Specifies the host name of the network location.
Port	Specifies the port number of the network location.
Protocol	Selects the protocol name of the network location.
Encoding	Selects the encoding option of the network location.
Test	Validates the connection to the network location.

To view the audit trail:

1. Access the security configuration settings.
2. Click the **Detail** selection for **Audit Trail**.
3. Click **Audit Trail View**.
 - **Audit Filter** lists the filters available for security log files.
 - **Audit Viewer** lists the viewing selections for security log files.

Archiving the Audit Trail

You can archive the audit trail by saving the security log file to the location specified in the configuration settings.

To archive the audit trail:

1. Access the security configuration settings.
2. Click the **Detail** selection for **Audit Trail**.
3. Click **Audit Trail Archive**.

The system displays an icon next to the name of the security log file to indicate the archive status of the log file.

Icon	Function
	The security log file is not stored.
	The security log file is stored.

4. Select the file and then click **Archive**.
The ultrasound system displays the disk space available on the system in megabytes and total quota limit in percentage.
5. To initiate the archive process after an incomplete archive, click **Restart**.
6. To remove the security log file from the ultrasound system, click **Delete**.

Appendix A EMC Publication

Electromagnetic Compatibility	2
Avoiding Electromagnetic Interference	2
Electrosurgical Units	3
Restrictions for Use.....	3
Immunity Level Test Results.....	3
Transducers, Accessories, and Cables	4
User-Accessible Cables.....	4

Electromagnetic Compatibility

The Electromagnetic compatibility (EMC) testing of this system was performed according to the international standard for EMC of medical devices (IEC 60601-1-2). This International Electrotechnical Commission (IEC) standard has been adopted in Europe as the European Norm standard (EN 60601-1-2). The ultrasound system from Siemens Healthineers meets the requirements of the standard and is suitable for a professional healthcare environment. Compliance with test limits does not guarantee that a particular installation will be free from electromagnetic interference.

The following is the definition for the essential performance of the ultrasound equipment:

- The image display is free of noise or distortion of the image or numerical value errors which may cause misdiagnosis.
- The image display is free from incorrect safety-related indications, for example, MI or TI.
- The device is free from unintended or excessive ultrasound output, transducer surface temperature, or uncontrolled motion of transducer assemblies intended for intra-corporeal (endocavity) use.

Avoiding Electromagnetic Interference

⚠ WARNING: Operating the ultrasound system in close proximity to other equipment can cause reciprocal interference. You should observe and ensure normal operation of the ultrasound system and other equipment.

Medical devices either generate (emit) or receive electromagnetic interference. The EMC standards describe tests for both emitted and received interference. Emission tests address interference generated by the device being tested. The ultrasound system from Siemens Healthineers meets the standard for radiated emissions.

Radio frequency (RF) is a form of electromagnetic energy tested by the EMC standards. Ultrasound systems and ultrasound transducers operate in the range of RF frequencies and are therefore susceptible to interference generated by other RF energy sources. Sources of RF interference include medical devices, information technology products, and radio/television transmission towers. To reduce interference from an RF energy source:

- Increase the separation (distance) between the ultrasound equipment and the interfering device (RF energy source).
- Connect the ultrasound equipment and the interfering device (RF energy source) to separate outlets on separate circuits, if possible.

In some environments, it is difficult to identify the source of radiated interference. To identify the source of the interference, determine the answers to these questions and then contact your service representative with the information you have collected.

- Is the interference intermittent or constant?
- Is the interference present with only one transducer or with several transducers?
- Do two different transducers operating at the same frequency have the same problem?
- Is the interference present if the ultrasound system is moved to a different location in the facility?
- Is the interference mode-dependent? That is, does the interference occur only in one mode, such as during Doppler?
- Can the EMC coupling path be attenuated? For example, does placing a transducer close to an ECG cable increase electromagnetic interference? Does moving the cable or other medical equipment away from the transducer result in reduced electromagnetic interference?

Electrosurgical Units

Electrosurgical units (ESUs) and other devices intentionally introduce radio frequency electromagnetic fields or currents into patients. Because ultrasound imaging frequencies are coincidentally in this radio frequency range, ultrasound transducer circuits are susceptible to radio frequency interference. When an ESU is in use, expect the ultrasound image to experience severe interference that may make the image unusable. This interference stops as soon as the ESU is powered off.

Restrictions for Use

If RF interference is detected, then the physician must determine if an artifact caused by the interference will negatively impact image quality and the subsequent diagnosis.

Immunity Level Test Results

Immunity is defined in the standard as the ability of a system to perform without degradation in the presence of electromagnetic disturbance. The EMC standards require manufacturers of patient-coupled equipment to specify immunity levels for their systems. The standards recognize that ultrasound equipment is designed to receive and amplify low-level signals in the same bandwidth as the interference. Therefore, it is reasonable to expect image noise during electromagnetic disturbance.

Electromagnetic Compatibility (EMC) Note: Operating the ultrasound imaging system in close proximity to sources of strong electromagnetic fields, such as radio transmitter stations or similar installations, as well as portable and mobile RF (radiated frequency) may lead to temporary degradation or interference visible on the monitor screen. A lightening of image background may be noticed while visualizing hypoechoic structures, or color spectral interference, or jitter, or horizontal lines in the image screen may occur. However, the device has been designed and tested to withstand such interference and will not be permanently damaged. Wait until the interference disappears or restart the ultrasound system.

When comparing immunity levels of different ultrasound systems, recognize that although the EMC standard does define the test methodology, it does not specify the criteria to assess degradation. Degradation assessment may vary by manufacturer.

A qualitative assessment of degradation in image quality is subjective. Noting when the first sign of an artifact is seen in the image minimizes the issue of subjectivity and also provides for stringent test results.

The results of emissions testing and immunity testing are provided in the accompanying tables. The standards used in testing are also provided in the tables. Testing was performed on a typical ultrasound system configuration. The values in the compliance level column of the tables specific to transducers is the level at which the RF interference is no longer discernable from the ambient background noise.

Transducers, Accessories, and Cables

⚠ WARNING: Accessory equipment connected to the analog and digital interfaces must be certified according to the respective EN or IEC standards (for example, EN 60950 or IEC 60950 for data processing equipment and EN 60601-1 or IEC 60601-1 for medical equipment). Anyone who connects additional equipment to any of the signal input or signal output ports configures a medical system and is therefore responsible that the system complies with the requirements of the system standards EN 60601-1 or IEC 60601-1. Siemens Healthineers can only guarantee the performance and safety of the devices listed in the *Instructions for Use*. If in doubt, consult the Siemens Healthineers service department or your local Siemens Healthineers representative.

⚠ WARNING: You must only use the transducers, accessories, cables, and replacement parts for internal components specified by Siemens Healthineers to reduce the risk of increased RF (radio frequency) emissions or decreased immunity of the ultrasound system.

User-Accessible Cables

See also: For information about cables, refer to Appendix A in the Instructions for Use.

Table 1. Guidance and manufacturer's declaration – electromagnetic emissions for the ultrasound system

The ultrasound system is intended for use in the electromagnetic environment specified below. The customer or the user of the ultrasound system should assure that it is used in such an environment.		
Emissions Test	Compliance	Electromagnetic Environment - Guidance
RF emissions CISPR 11	Group 1	The ultrasound system uses RF energy only for its internal function. Therefore, its RF emissions are very low and are not likely to cause any interference in nearby electronic equipment.
RF emissions CISPR 11	Class B Note: The system is a Class A device when the barcode reader, DVR, or 9VC2, 10BCC3, or 9VE4 transducer is in use.	The ultrasound system is suitable for use in all establishments other than domestic and those directly connected to the public low-voltage power supply network that supplies buildings used for domestic purposes.
Harmonic emissions IEC 61000-3-2	Class A	
Voltage fluctuations/flicker emissions IEC 61000-3-3	Complies	

Table 2. Guidance and manufacturer's declaration – electromagnetic immunity for the ultrasound system independent of transducer type

The ultrasound system is intended for use in the electromagnetic environment specified below. The customer or the user of the ultrasound system should assure that it is used in such an environment.			
Immunity Test	IEC 60601-1-2 Test Level	Compliance Level	Electromagnetic Environment - Guidance
Electrostatic discharge (ESD) IEC 61000-4-2	±8 kV contact ±15 kV air	±8 kV contact ±15 kV air	Floors should be wood, concrete or ceramic tile. If floors are covered with synthetic material, the relative humidity should be at least 30%.
Electrical fast transient/burst IEC 61000-4-4	±2 kV for power supply lines ±1 kV for input/output lines	±2 kV for power supply lines ±1 kV for input/output lines	Mains power quality should be that of a typical commercial or hospital environment.
Surge IEC 61000-4-5	±1 kV line to line ±2 kV line to ground	±1 kV differential mode ±2 kV common mode	Mains power quality should be that of a typical commercial or hospital environment.
Voltage dips, short interruptions and voltage variations on power supply input lines IEC 61000-4-11	0% U_T (100% dip in U_T) for 0.5 cycle	System may shut down and the mains circuit breaker could be switched to the off position (O = OFF).	Ensure the mains circuit breaker is in the on position (I = ON) and then use the power on sequence to reboot the system. Mains power quality should be that of a typical commercial or hospital environment. If the user of the ultrasound system requires continued operation during power mains interruptions, it is recommended that the ultrasound system be powered from an uninterruptible power supply or a battery.
	0% U_T (100% dip in U_T) for 1 cycle	0% U_T (100% dip in U_T) for 1 cycle	Mains power quality should be that of a typical commercial or hospital environment. If the user of the ultrasound system requires continued operation during power mains interruptions, it is recommended that the ultrasound system be powered from an uninterruptible power supply or a battery.
	70% U_T (30% dip in U_T) for 25 cycles	70% U_T (30% dip in U_T) for 25 cycles	
	0% U_T (100% dip in U_T) for 250 cycles	System goes to normal power down condition.	For extended mains power outages without the use of an uninterruptible power supply or a battery, the system will shut down normally. Use the power on sequence to reboot the system.
Power frequency (50/60 Hz) magnetic field IEC 61000-4-8	30 A/m	30 A/m	Power frequency magnetic fields should be at levels characteristic of a typical location in a typical commercial or hospital environment.
Note: U_T is the a.c. mains voltage prior to application of the test level.			

Table 3. Guidance and manufacturer's declaration – electromagnetic immunity for the ultrasound system

The ultrasound system is intended for use in the electromagnetic environment specified below. The customer or the user of the ultrasound system should assure that the system is used in such an environment.			
Immunity Test	IEC 60601-1-2 Test Level	Compliance Level	Electromagnetic Environment – Guidance
Conducted RF IEC 61000-4-6	3 Vrms 150 kHz to 80 MHz outside ISM bands	3 Vrms	Portable and mobile RF communications equipment should be used no closer to any part of the ultrasound system, including cables, than the recommended separation distance calculated from the equation applicable to the frequency of the transmitter. Recommended separation distance (<i>d</i>) $d = 1.2 \sqrt{P}$ $d = 2.0 \sqrt{P}$ where <i>P</i> is the maximum output power rating of the transmitter in watts (W) according to the transmitter manufacturer and <i>d</i> is the recommended separation distance in meters (m). Field strengths from fixed RF transmitters, as determined by an electromagnetic site survey, ^a should be less than the compliance level in each frequency range. ^b Interference may occur in the vicinity of equipment marked with the following symbol: 
Radiated RF IEC 61000-4-3	6 Vrms 150 kHz to 80 MHz within ISM bands	6 Vrms	
	3 V/m 80 MHz to 2.7 GHz	3 V/m	
Note 1: At 80 MHz and 800 MHz, the higher frequency range applies.			
Note 2: These guidelines may not apply in all situations. Electromagnetic propagation is affected by absorption and reflection from structures, objects, and people.			
Note 3: The ISM (Industrial, Scientific, and Medical) bands between 150 kHz and 80 MHz are 6.765 MHz to 6.795 MHz; 13.553 MHz to 13.567 MHz; 26.957 MHz to 27.283 MHz; and 40.66 MHz to 40.70 MHz.			
^a Field strengths from fixed transmitters, such as base stations for radio (cellular/cordless) telephones and land mobile radios, amateur radio, AM and FM radio broadcast, and TV broadcast cannot be predicted theoretically with accuracy. To assess the electromagnetic environment due to fixed RF transmitters, an electromagnetic site survey should be considered. If the measured field strength in the location in which the ultrasound system is used exceeds the applicable RF compliance level above, the ultrasound system should be observed to verify normal operation. If abnormal performance is observed, additional measures may be necessary, such as changing transducers or transducer operating frequency, or re-orienting or relocating the ultrasound system.			
^b Over the frequency range 150 kHz to 80 MHz, field strengths should be less than the compliance level 3 V/m.			

Table 4. Recommended separation distances between portable and mobile RF communications equipment and the ultrasound system

The ultrasound system is intended for use in an electromagnetic environment in which radiated RF disturbances are controlled. The customer or the user of the ultrasound system can help prevent electromagnetic interference by maintaining a minimum distance between portable and mobile RF communications equipment (transmitters) and the ultrasound system as recommended below, according to the maximum output power of the communications equipment.		
Rated maximum output power of transmitter W (watts)	Separation distance according to frequency of transmitter m (meters)	
	150 kHz to 80 MHz	80 MHz to 2.7 GHz
	$d = 1.2 \sqrt{P}$	$d = 2.0 \sqrt{P}$
0.01	0.12	0.20
0.1	0.38	0.63
1	1.2	2.0
10	3.8	6.3
100	12	20
For transmitters rated at a maximum output power not listed above, the recommended separation distance d in meters (m) can be estimated using the equation applicable to the frequency of the transmitter, where P is the maximum output power rating of the transmitter in watts (W) according to the transmitter manufacturer.		
Note 1: At 80 MHz and 800 MHz, the separation distance for the higher frequency range applies.		
Note 2: These guidelines may not apply in all situations. Electromagnetic propagation is affected by absorption and reflection from structures, objects, and people.		

Table 5. Guidance and manufacturer's declaration - electromagnetic immunity

The ultrasound system is intended for use in an electromagnetic environment in which radiated RF disturbances are controlled. Portable RF communications equipment should be used no closer than 30 cm (12 inches) to any part of the ultrasound system. Otherwise, degradation of the performance of this equipment could result.					
Immunity Test	Band ^a	Service ^a	Modulation	IEC 60601 test level	Compliance level
Proximity fields from RF wireless communications equipment IEC 61000-4-3	380 - 390 MHz	TETRA 400	Pulse modulation 18 Hz	27 V/m	27 V/m
	430 - 470 MHz	GMRS 460 FRS 460	FM ± 5 kHz deviation 1 kHz sine	28 V/m	28 V/m
	704 - 787 MHz	LTE Band 13, 17	Pulse modulation 217 Hz	9 V/m	9 V/m
	800 - 960 MHz	GSM 800/900 TETRA 800 iDEN 820 CDMA 850 LTE Band 5	Pulse modulation 18 Hz	28 V/m	28 V/m
	1700 - 1990 MHz	GSM 1800 CDMA 1900 GSM 1900 DECT LTE Band 1, 3, 4, 25 UMTS	Pulse modulation 217 Hz	28 V/m	28 V/m
	2400 - 2570 MHz	Bluetooth WLAN 802.11 b/g/n RFID 2450 LTE Band 7	Pulse modulation 217 Hz	28 V/m	28 V/m
	5100 - 5800 MHz	WLAN 802.11 a/n	Pulse modulation 217 Hz	9 V/m	9 V/m
Note: These guidelines may not apply in all situations. Electromagnetic propagation is affected by absorption and reflection from structures, objects, and people.					
^a For some services, only the uplink frequencies are included.					

Appendix B Measurements and Calculations Reference

Measurement Methods	2
2D-Mode Measurement Methods	2
M-Mode Measurement Methods	3
Doppler Measurement Methods	4
Exam-Specific Measurement Methods	6
2D-Mode Obstetric Measurement Methods	6
Measurement Formulas	7

Measurement Methods

The following general measurements apply to all exam types.

2D-Mode Measurement Methods

Measurement Method	Measurement Results	Displayed Results
Distance	Distance	D cm
4 Point Dist	Distance	D1, D2, D3 cm
5 Point Dist	Distance	D1, D2, D3, D4 cm
6 Point Dist	Distance	D1, D2, D3, D4, D5 cm
2 Distance Avg	Distance	D1, D2 cm Avg cm
3 Distance Avg	Distance	D1, D2, D3 cm Avg cm
Open Trace	Distance	D cm
Open Spline	Distance	D cm
Ellipse	Distance	D1, D2 cm
	Circumference	C cm
	Area	A cm ²
Single Angle Cardiac Axis	Angle	Angle °
Hip (α)	Angle	α °
Hip (α+β)	Angle	α °
Hip (β+α)		β °
Cross Angle		
1 Ellipse	Distance	D1, D2 cm
	Circumference	C cm
	Area	A cm ²
	Volume	V cm ³
Vol 1E + 1D	Distance	D, D1, D2 cm
	Circumference	C cm
	Area	A cm ²
	Volume	V cm ³
Vol 1 Dist	Distance	D cm
	Volume	V cm ³
Vol 2 Dist	Distance	D1, D2 cm
	Volume	V cm ³
Vol 3 Dist	Distance	D1, D2, D3 cm
	Volume	V cm ³

Measurement Method	Measurement Results	Displayed Results
Disk	Disk Volume	V cm³
Disk Trace		
Dist Ratio	Distance	A, B cm
Volume Ratio	Ratio	A/B
Area Ratio		B/A
d:D Ratio		(A-B)/A
%FHC	Distance	D cm
	Distance	d cm
	Ratio	%FHC %
eSieCalcs	Distance	D cm
	Circumference	C cm
	Area	A cm²
Trace	Circumference	C cm
Spline	Area	A cm²
D-%Stenosis	Distance	D1, D2 cm
	Percentage stenosis	%Steno %
A-%Stenosis	Area	A1, A2 cm²
ET-%Stenosis	Percentage stenosis	%Steno %
Auto Stenosis	Area	A1, A2 cm²
	Circumference	C1, C2 cm
	Percentage stenosis	%Steno %
O-Marker	Distance	D cm

M-Mode Measurement Methods

Measurement Method	Measurement Results	Displayed Results
Distance	Distance	D cm
4 Point Dist	Distance	D1, D2, D3 cm
5 Point Dist	Distance	D1, D2, D3, D4 cm
Slope	Distance	D cm
	Slope	Slope cm/s
	Time	T ms
HR	Heart rate	HR bpm
	Time	T ms
Time	Time	T ms
Dist Ratio	Distance	D cm
	Ratio	A/B
		B/A
		(A-B)/A
Time Ratio	Time	T s
	Ratio	A/B
		B/A
		(A-B)/A

Doppler Measurement Methods

Measurement Method	Measurement Results	Displayed Results
PI Auto	Pulsatility index	PI
	Resistive index	RI
	Systole/diastole	S/D
	Peak systole	PS cm/s
	End diastole	ED cm/s
	Time-averaged maximum	TAmx
	Time-averaged mean	TAmx
PI Manual	Pulsatility index	PI
	Resistive index	RI
	Systole/diastole	S/D
	Peak systole	PS cm/s
	End diastole	ED cm/s
	Time-averaged maximum	TAmx cm/s
Velocity	Velocity	V cm/s
	Frequency	Freq kHz
	Pressure gradient	PG mmHg
RI	Peak systole	PS cm/s
	End diastole	ED cm/s
	Systole/diastole	S/D
	Resistive index	RI
A-Flow	Distance	D1, D2 cm
	Circumference	C cm
	Area	A cm²
	Time-averaged mean	TAMean cm/s
	Flow volume	FV L/min
D-Flow	Distance	D1 cm
	Area	A cm²
	Time-averaged mean	TAMean cm/s
	Flow volume	FV L/min
Velocity Ratio	Velocity	V1, V2 cm/s
	Frequency	Freq1, Freq2 kHz
	Ratio	V1/V2
	Ratio	V2/V1
	Ratio	(V1-V2)/V1

Measurement Method	Measurement Results	Displayed Results
Acceleration	Slope	Acceleration m/s²
	Time	T ms
	Peak velocity	Vmx cm/s
Time	Time	T ms
VTI	Velocity-time integral	VTI cm
	Maximum velocity	Vmax m/s
	Mean velocity	Vmean m/s
	Maximum pressure gradient	PGmax mmHg
	Mean pressure gradient	PGmean mmHg
	Actual time	AT ms
	Estimated time	ET ms
VTI Trace	Velocity-time integral	VTI cm
	Maximum velocity	Vmax m/s
	Mean velocity	Vmean m/s
	Maximum pressure gradient	PGmax mmHg
	Mean pressure gradient	PGmean mmHg
	Actual time	AT ms
	Estimated time	ET ms
HR	Heart rate	HR bpm
	Time	T ms
PHT Guide	Pressure half time	PHT ms
	Maximum velocity	Vmax cm/s
	Velocity	V1, V2 cm/s
	Deceleration time	DT ms
	Deceleration slope	DS m/s²
PHT	Pressure half time	PHT ms
	Maximum velocity	cm/s
	Velocity	V1, V2 cm/s
	Deceleration time	DT ms
	Deceleration slope	DS m/s²
dP/dt	Delta time	dt ms
	Delta pressure over delta time	dP/dt mmHg/s
Time Ratio	Time	T s
	Ratio	A/B
		B/A
		(A-B)/A
E/A Vel + DT	Pressure half time	PHT ms
	Deceleration time	DT ms
	Deceleration slope	DS m/s²
	E point	E Pt cm/s
	A point	A Pt cm/s

Exam-Specific Measurement Methods

2D-Mode Obstetric Measurement Methods

The following measurement methods are available only during the following exam types:
obstetric, early obstetric.

Measurement Method	Measurement Results	Displayed Results
AFI	Distance	AFI mm
		Q1, Q2, Q3, Q4 mm
AP Depth	Distance	D cm

Measurement Formulas

2D-Mode Measurement Formulas

Measurement Result	Formula	Abbreviation
Area	$A = \pi \times D1 / 2 \times D2 / 2 / 100 \text{ cm}^2$ $C = \pi (D1/2 + D2/2) \text{ mm}$	A = area C = circumference
Percentage stenosis, distance	$\%Steno = (D1 - D2) / D1 \times 100$	%Steno = percentage stenosis D1 = distance of long axis D2 = distance of short axis
Percentage stenosis, area	$\%Steno = (A1 - A2) / A1 \times 100$	%Steno = percentage stenosis A1 = large area A2 = small area
Volume, one distance	$V = \pi \times D^3 / 6 / 1000 \text{ cm}^3$	V = volume D = distance
Volume, two distances	$V = \pi \times D1 \times D2^2 / 6 / 1000 \text{ cm}^3$	V = volume D1 = length D2 = distance
Volume, three distances	$V = \pi \times D1 \times D2 \times D3 / 6 / 1000 \text{ cm}^3$	V = volume D1 = length D2 = depth D3 = width
Volume, one ellipse and one distance	$C = \pi (D1/2 + D2/2) \text{ mm}$ $A = \pi \times D1\frac{1}{2} \times D2\frac{1}{2} / 100 \text{ cm}$ $V = \pi \times D1 \times D2 \times D/6 / 1000 \text{ cm}^3$	V = volume A = area C = circumference D1 = distance D2 = distance
Volume, one ellipse	$C = \pi (D1 / 2 + D2 / 2) \text{ mm}$ $A = \pi \times D1\frac{1}{2} \times D2\frac{1}{2} / 100 \text{ cm}$ $V = 8A^2 / (3 \times \pi (\text{longer value between } D1 \text{ \& } D2)) \times 10$	V = volume A = area C = circumference D1 = distance D2 = distance D = distance
Volume, disk	$V = \pi \times L / 20 \times \text{Sigma } (Ri^2 \text{ mm}) / 1000 \text{ mL}$	V = volume
Average, two distances	$Avg = (D1 + D2) / 2$	D1 = distance D2 = distance Avg: average of distances
Average, three distances	$Avg = (D1 + D2 + D3) / 3$	D1 = distance D2 = distance D3 = distance Avg: average of distances
Angle, cross angle	$\alpha + \beta = 180^\circ$ $\alpha < 180^\circ, \beta < 180^\circ$	α = angle β = angle

M-Mode Measurement Formulas

Measurement Result	Formula	Abbreviation
Slope	$S = D / T$	S = slope D = distance T = time

Doppler Measurement Formulas

Measurement Result	Formula	Abbreviation
Resistive index	$RI = (PS - ED) / PS$	RI = resistive index PS = peak systole ED = end diastole
Pulsatility index (end diastole)	$PI (ED) = (PS - ED) / TAmx$	PI (ED) = pulsatility index (end diastole) PS = peak systole ED = end diastole TAmx = time averaged maximum velocity
Pulsatility index (minimum diastole)	$PI (MD) = (PS - MD) / TAmx$	PI (MD) = pulsatility index (minimum diastole) PS = peak systole MD = minimum diastole TAmx = time averaged maximum velocity
Pulsatility index (transcranial Doppler)	$PI (TCD) = 3 \times (PS - ED) / \{PS + (2 \times ED)\}$	PI (TCD) = pulsatility index (transcranial Doppler) PS = peak systole ED = end diastole
Mean Transcranial Doppler	$Mean\ TCD = \{PS + (2 \times ED)\} / 3$	Mean TCD = mean transcranial Doppler PS = peak systole ED = end diastole
Acceleration	$Accel = (Vmx - Vmn) \text{ cm} / T \text{ ms} \times 1000 \text{ sec} / 100 \text{ m}$	V = volume T = time Accel = acceleration
Flow volume	$FV = TAmn \times A \times 60 / 1000$ $(A = \pi / 4 \times D^2 / 100)$	FV = flow volume TAmn = time averaged mean velocity A = area D = distance
Flow volume	$FV = TAmn \times A \times 60 / 1000$ $(A = \pi / 4 \times D1 \times D2 / 100)$	FV = flow volume TAmn = time averaged mean velocity A = area D1, D2 = distance
Velocity	$P = 4V^2 \text{ m}$	P = peak gradient V = volume
Velocity time integral	$VTI = \text{Velocity} \times \text{Time}$	VTI = velocity time integral
Heart rate	$\{60 \text{ sec} / 1 \text{ heart beat cycle time} \times 1000 \text{ sec}\}$	-

See also: For information about exam specific reports, refer to Appendix C, Appendix D, and Appendix E in this manual.

Appendix C Obstetrical Reference

Obstetric Measurement and Calculation Labels.....	2
2D-Mode Measurement Labels for Estimating Gestational Age	2
2D-Mode Measurement Labels Not Related to Estimating Gestational Age	3
2D-Mode Calculation Labels.....	6
Doppler Calculation Labels.....	7
M-Mode Measurement Labels	8
Doppler Measurement Labels.....	8
M-Mode Calculation Labels	9
Sections in the Patient Report.....	10
Growth Analysis Graph Measurements and Calculations.....	11

Obstetric Measurement and Calculation Labels

2D-Mode Measurement Labels for Estimating Gestational Age

Note: An estimated fetal age can be calculated when you measure both **APTD** and **TTD**.

Measurement Label	Description	Measurement Method
Mean Sac Diam	Mean Gestational Sac Diameter	3 Distance Avg
GSD	Gestational Sac, Maximum Length	Distance
CRL	Crown Rump Length	Distance
BPD	Biparietal Diameter	Distance
HC	Head Circumference. Also performs the following measurements and calculations: ▪ Cephalic Index (CI) ▪ Biparietal Diameter (BPD), if specified in the configuration settings ▪ Occipital Frontal Diameter (OFD), if specified in the configuration settings	Ellipse
AC	Abdominal Circumference. Also performs the following measurements and calculations if specified in the configuration settings: ▪ Antero-posterior Abdominal Diameter (APAD) ▪ Transverse Abdominal Diameter (TAD) Note: When you measure abdominal circumference, the system assigns the larger measurement to transverse abdominal diameter and the smaller measurement to antero-posterior abdominal diameter.	Ellipse
FL	Femur Length	Distance
OFD	Occipital Frontal Diameter	Distance
TAD	Transverse Abdominal Diameter	Distance
APAD	Antero-posterior Abdominal Diameter	Distance
HL	Humerus Length	Distance
Binoc D	Binocular Distance (Outer Orbital Diameter)	Distance
Clavicle	Clavicle Length	Distance
FTA	Fetal Trunk Area	Ellipse
APTD	Antero-posterior Trunk Diameter	Distance
TTD	Transverse Thoracic Diameter	Distance
Tibia	Tibia Length	Distance
Ulna	Ulna Length	Distance
Foot	Foot Length	Distance

2D-Mode Measurement Labels Not Related to Estimating Gestational Age

Measurement Label	Description	Measurement Method
Yolk Sac	Yolk Sac Diameter	Distance
Intra Trans	Intracranial Translucency	Distance
Nuchal Translucency	Nuchal Translucency	Distance
Nuchal Thickness	Nuchal Fold Thickness	Distance
TC	Thoracic Circumference	Ellipse
AFI	Amniotic Fluid Index in up to four quadrants	Distance
HW	Hemispheric Width (Width of Hemisphere)	Distance
LV-A	Anterior Horn Lateral Ventricular Width	Distance
LV-P	Posterior Horn Lateral Ventricular Width	Distance
Lat Vent	Lateral Ventricle Width	Distance
Cist Magna	Cisterna Magna	Distance
Cervix (TA)	Cervix Length, transabdominal	Distance
Cervix (EV)	Cervical Length, endovaginal	Distance
Umb V Diam	Umbilical Vein Diameter	Distance
APTD	Anteroposterior Trunk Diameter	Distance
TTD	Transverse Trunk Diameter	Distance
Facial Angle	Facial Angle	Single Angle
TA	Thoracic Area	Ellipse
MAD	Mean Abdominal Diameter	Distance
Fibula	Fibula Length	Distance
Radius	Radius Length	Distance
Kidney Length	Kidney Length	Distance
Kidney Width	Kidney Width	Distance
Kidney AP	Kidney Anteroposterior	Distance
Ureter	Ureter Length	Distance
PreVoid Bladder Length	Pre-Void Bladder Length	Distance
PreVoid Bladder AP	Pre-Void Bladder Anteroposterior	Distance
PreVoid Bladder Width	Pre-Void Bladder Width	Distance
PostVoid Bladder Length	Post-Void Bladder Length	Distance
PostVoid Bladder AP	Post-Void Bladder Anteroposterior	Distance
PostVoid Bladder Width	Post-Void Bladder Width	Distance
Fibroid 1	Fibroid 1, numbered 1 through 5	Distance
Cyst 1	Cyst 1, numbered 1 through 5	Distance
Mass 1	Mass 1, numbered 1 through 5	Distance

Measurement Label	Description	Measurement Method
CTR Area	Cardiothoracic Area Ratio	Area Ratio
CSP	Cavum Septum Pellucidum	Distance
IOD	Interocular Distance (Inner Ocular Distance)	Distance
Nasal Bone	Nasal Bone	Distance
Cerebellum	Cerebellum (Trans Cerebellar Diameter)	Distance
Aorta	Aorta Diameter	Distance
Umb A	Umbilical Artery	Distance
Cardiac Axis	Cardiac Axis	Cardiac Axis
Amniotic Pocket	Amniotic fluid pocket	Distance
Q1	Amniotic Fluid Index Q1, numbered Q1 through Q4	Distance
Uterus Width (TA)	Uterus Width, transabdominal	Distance
Uterus Length (TA)	Uterus Length, transabdominal	Distance
Uterus AP (TA)	Uterus Anteroposterior Diameter, transabdominal	Distance
Uterus Width (EV)	Uterus Width, endovaginal	Distance
Uterus Length (EV)	Uterus Length, endovaginal	Distance
Uterus AP (EV)	Uterus Anteroposterior Diameter, endovaginal	Distance
Ovary Width (TA)	Ovary Width, transabdominal	Distance
Ovary Length (TA)	Ovary Length, transabdominal	Distance
Ovary AP (TA)	Ovary Anteroposterior Diameter, transabdominal	Distance
Ovary Width (EV)	Ovarian Transverse Diameter (Width), endovaginal	Distance
Ovary Length (EV)	Ovary Length, endovaginal	Distance
Ovary AP (EV)	Ovary Anteroposterior Diameter, endovaginal	Distance
Endo Thick Width (TA)	Endometrium Thickness Width, transabdominal	Distance
Endo Thick Length (TA)	Endometrium Thickness Length, transabdominal	Distance
Endo Thick AP (TA)	Endometrium Thickness Anteroposterior Diameter, transabdominal	Distance
Endo Thick Width (EV)	Endometrium Thickness Width, endovaginal	Distance
Endo Thick Length (EV)	Endometrium Thickness Length, endovaginal	Distance
Endo Thick AP (EV)	Endometrium Thickness Anteroposterior Diameter, endovaginal	Distance

Measurement Label	Description	Measurement Method
AoD Arch	Aortic Arch Diameter	Distance
LVOT diam	Left Ventricular Outflow Tract	Distance
RVOT diam	Right Ventricular Outflow Tract	Distance
Asc AoD	Ascending Aorta Diameter	Distance
Desc AoD	Descending Aorta Diameter	Distance
PA	Pulmonary Artery	Distance
IVS d	Interventricular Septal Thickness in diastole	Distance
IVS s	Interventricular Septal Thickness in systole	Distance
LVID d	Left Ventricle Internal Diameter in diastole	Distance
LVID s	Left Ventricle Internal Diameter in systole	Distance
LVPW d	Left Ventricle Posterior Wall Thickness in diastole	Distance
LVPW s	Left Ventricle Posterior Wall Thickness in systole	Distance
RVD d	Right Ventricle Diameter in diastole	Distance
Total Cardiac Dimension	Total Cardiac Dimension	Distance

2D-Mode Calculation Labels

Calculation Label	Description	Required Measurements
EFW1 EFW2	Estimated Fetal Weight	Required measurements vary depending on the selected reference authors
PreVoid Bladder Vol	Pre-Void Bladder Volume	PreVoid Bladder Length PreVoid Bladder AP PreVoid Bladder Width
PostVoid Bladder Vol	Post-Void Bladder Volume	PostVoid Bladder Length PostVoid Bladder AP PostVoid Bladder Width
Micturated Volume	Micturated Volume	PreVoid Bladder Vol PostVoid Bladder Vol
Ovary Volume (TA)	Ovary, transabdominal (Transabdominal)	Ovary Length (TA) Ovary Width(TA) Ovary AP (TA)
Ovary Volume (EV)	Transvaginal Ovarian Volume (Endovaginal)	Ovary Length (EV) Ovary Width(EV) Ovary AP (TA)
AFI	Amniotic Fluid Index	Q1 Q2 Q3 Q4
Uterus Volume (TA)	Uterus Volume (Transabdominal)	Uterus (TA) Length Uterus (TA) Width Uterus (TA) AP
Uterus Volume (EV)	Uterus Volume (Endovaginal)	Uterus (EV) Length Uterus (EV) Width Uterus (EV) AP
Kidney Volume (Fetus)	Kidney Volume (Fetus)	Kidney Length Kidney Width Kidney AP
APTD X TTD	Anteroposterior Trunk Diameter multiplied by Transverse Trunk Diameter	APTD TTD
BPDa	Biparietal Diameter corrected by a formula to adjust the head shape	BPD OFD
Cerebellum/ AC	Growth Analysis from Transverse Cerebellar Diameter over Abdominal Circumference	Cerebellum AC

Calculation Label	Description	Required Measurements
FL/AC	Ratio of Femur Length to Abdominal Circumference	FL AC
FL/BPD	Ratio of Femur Length to Biparietal Diameter	FL BPD
HC/AC	Ratio of Head Circumference to Abdominal Circumference	HC AC
Lat Vent/HW	Ratio of Lateral Ventricular Width to Hemispheric Width	LVW HW
CI	Cephalic Index short axis / long axis * 100	BPD OFD
Cubed / Gibson / Teichholz	---	
EDV	Left Ventricular End-Diastolic Volume	LVID d
EF	Ejection Fraction is the ratio of the stroke volume to the end-diastolic volume	SV EDV
ESV	Left Ventricular End-Systolic Volume	LVID s
FS	Fractional Shortening is the percentage of shortening of the Left ventricular dimension	LVID d LVID s

Doppler Calculation Labels

Calculation Label	Description	Required Measurements
Ductus V PLI	Ductus Venous Preload Index	Ductus V S Vmax Ductus V A Vmax
Ductus V PVIV	Ductus Venous Peak velocity index for veins	Ductus V S Vmax Ductus V A Vmax Ductus V D Vmax
Ductus V S/a	Ductus Venous S point / a point	Ductus V S Vmax Ductus V A Vmax
LV MPI	Left Ventricle Myocardial Performance Index	LV ET LV C-O dur
RV MPI	Right Ventricle Myocardial Performance Index	RV ET RV C-O dur

M-Mode Measurement Labels

Measurement Label	Description	Measurement Method
Fetal HR	Fetal Heart Rate, in beats per minute	HR
LVOT diam	Left Ventricular Outflow Tract	Distance
RVOT diam	Right Ventricular Outflow Tract	Distance
Asc AoD	Ascending Aorta Diameter	Distance
Desc AoD	Descending Aorta Diameter	Distance
AoD Arch	Aortic Arch Diameter	Distance
PA	Pulmonary Artery	Distance
IVS d	Interventricular Septal Thickness in diastole	Distance
IVS s	Interventricular Septal Thickness in systole	Distance
LVID d	Left Ventricle Internal Diameter in diastole	Distance
LVID s	Left Ventricle Internal Diameter in systole	Distance
LVPW d	Left Ventricle Posterior Wall Thickness in diastole	Distance
LVPW s	Left Ventricle Posterior Wall Thickness in systole	Distance
RVD d	Right Ventricle Diameter in diastole	Distance

Doppler Measurement Labels

Measurement Label	Description	Measurement Method
Fetal HR	Fetal Heart Rate in beats per minute	HR
Umb A	Umbilical Artery	PI Auto
Umb V A Vmax	Umbilical Vein A Vmax	Velocity
Umb V D Vmax	Umbilical Vein D Maximum Velocity	Velocity
Umb V S Vmax	Umbilical Vein S Maximum Velocity	Velocity
Uterine A	Uterine Artery	PI Auto
Ductus V All	Ductus Venous All	PI Auto
Ductus V S	Ductus V S Vmax	Velocity
Ductus V A	Ductus V A Vmax	Velocity
Ductus V D	Ductus V D Vmax	Velocity
Aorta	Aorta	PI Auto
IVC PLI	Preload Index	Velocity Ratio V1 = S, V2=a
Ao Desc s	Aortic Descending Systole	Velocity

Measurement Label	Description	Measurement Method
LVOT Vmax	Left Ventricular Outflow Tract Maximum Velocity	Velocity
RVOT Vmax	Right Ventricular Outflow Tract Maximum Velocity	Velocity
LV ET	Left Ventricular Ejection Time	Time
RV ET	Right Ventricular Ejection Time	Time
LV C-O dur	Left Ventricular Close-Open Duration	Time
RV C-O dur	Right Ventricular Close-Open Duration	Time
Asc Ao Vmax	Ascending Aorta Maximum Velocity	Velocity
Desc Ao Vmax	Descending Aorta Maximum Velocity	Velocity
DA Vmax	Ductus Arteriosus Maximum Velocity	Velocity
PA	Pulmonary Artery	Velocity
Ao Arch Vmax	Aortic Arch Velocity Maximum	Velocity
Ovarian A	Ovarian Artery	PI Auto
Uterine A	Uterine Artery	PI Auto
MCA	Middle Cerebral Artery	PI Auto

M-Mode Calculation Labels

Calculation Label	Description	Required Measurements
Cubed / Gibson / Teichholz	---	
EDV	Left Ventricular End-Diastolic Volume	LVID d
EF	Ejection Fraction is the ratio of the stroke volume to the end-diastolic volume	SV EDV
ESV	Left Ventricular End-Systolic Volume	LVID s
FS	Fractional Shortening is the percentage of shortening of the left ventricular dimension	LVID d LVID s
SV	Stroke Volume is the volume of blood ejected from a ventricle during one cardiac cycle or phase of ventricular systole	EDV ESV

Sections in the Patient Report

Section	Description
Exam Specific Information	Gestational age, date of confinement, and fetal weight estimations
2D Measurements	Measurement labels and results for 2D-mode measurements
2D Ratio	2D-mode measurement ratios specified in the configuration settings.
M Measurements	Measurement labels and results for M-mode measurements
D Measurements	Measurement labels and results for Doppler measurements
OB (2D), (D)	(Available only for the obstetric exam type) Measurement labels and result only for the obstetric exam in 2D or Doppler mode.
Fetal Heart (2D), (M), (D)	(Available only for the advanced obstetric exam types) Measurement labels and result for fetal heart motion in 2D, M or Doppler mode.
Assessment	Descriptions for exam-related structures. You can select system defined or user-defined descriptions. ▪ Fetal Heart Motion ▪ 4Ch heart ▪ RVOT ▪ LVOT ▪ Ao Arch ▪ Ductal Arch ▪ Fetal Number ▪ Fetal Position ▪ Placenta Location ▪ Previa ▪ Placenta Grade ▪ AF Volume ▪ Placental Cord ▪ Facial Profile ▪ Lips/Palate ▪ Lateral Ventricles ▪ Choroid ▪ Posterior Fossa ▪ Cerebellum ▪ Nuchal Fold ▪ C-Spine ▪ T-Spine ▪ L-Spine ▪ Sacrum ▪ 4 Extremities ▪ L Ankle/Foot ▪ Diaphragm ▪ Lungs ▪ Stomach ▪ Cord Insertion ▪ 3-V Cord ▪ Abdominal Wall ▪ R Kidney ▪ L Kidney ▪ Bladder ▪ Bowel Pattern ▪ Genitalia ▪ Cardiac Rhythm ▪ R Ankle/Foot

Section	Description
Biophysical Profile	Fetal parameters for creating a biophysical profile. ▪ Fetal Breathing Movements ▪ Fetal Tone ▪ Fetal Movements ▪ Amniotic Fluid Volume ▪ Total Score
Growth Graph	(Available only for the early obstetric and obstetric exam types) Growth analysis graphs indicating predicted fetal growth patterns
Quick Grow Graph	(Available only for the advanced obstetric exam type) Measurement and calculation results and quick growth graphs
USGA Data Graph	(Available only for the advanced obstetric exam type) Graphs with an estimated gestational age
Images	Images included from the thumbnail panel
Comments	Any combination of predefined comments and user-entered comments

Growth Analysis Graph Measurements and Calculations

References for growth analysis often differ from the references for estimated gestational age for the same label.

Note: References for growth curves typically use 5 or 10 percent for the lower limit and 90 or 95 percent for the upper limit.

The JSUM reference uses 5 percent and 95 percent for the lower and upper limits, respectively, for the middle cerebral artery and umbilical artery growth curves. The JSUM, Tokyo, and Osaka references use a standard deviation of 1.5 from the mean for the lower and upper limits for all other growth curves.

Graph	Description	References
Mean Sac Diam	Mean Gestational Sac Diameter	Hellman, Rempen
CRL	Crown Rump Length	Hadlock, Robinson, Hansmann, Rempen, ASUM, JSUM, Osaka, Tokyo, INTERGROWTH-21st, Pexters
BPD	Biparietal Diameter	Hadlock, Merz, Lasser, Hansmann, Rempen, ASUM, Chitty (O-I), Chitty (O-O), JSUM, Osaka, Tokyo, Nicolaides, Verburg, CFEF, Grumbach (Twin), INTERGROWTH-21st
OFD	Occipital Frontal Diameter	Hansmann, ASUM, Chitty, Nicolaides, INTERGROWTH-21st
HC	Head Circumference	Hadlock, Merz, Hansmann, ASUM, Chitty, Nicolaides, Verburg, CFEF, Deter, INTERGROWTH-21st
AC	Abdominal Circumference	Hadlock, Merz, Jeanty, ASUM, Chitty, JSUM, Nicolaides, Verburg, CFEF, Deter, Grumbach (Twin), INTERGROWTH-21st

Graph	Description	References
FL	Femur Length	Jeanty, Hadlock, Merz, Hansmann, ASUM, JSUM, Osaka, Tokyo, Nicolaides, Chitty, Verburg, CFEF, INTERGROWTH-21st, Warda
HL	Humerus Length	Jeanty, Merz, Hansmann, ASUM, Osaka, Chitty
Ulna	Ulna Length	Jeanty, Merz, Hansmann, Chitty
Tibia	Tibia Length	Jeanty, Merz, Hansmann, Chitty
Foot	Foot Length	Mercer, Chitty
TC	Thoracic Circumference	Chitkara
Clavicle	Clavicle Length	Yarkoni
FTA	Fetal Trunk Area	Osaka
AFI	Amniotic Fluid Index	Moore
Cerebellum	Cerebellum (Trans Cerebellar Diameter)	CFEF, Nicolaides, Verburg
Cist Magna	Cisterna Magna (Posterior Fossa)	Nicolaides
EFW1, EFW2	Estimated Fetal Weight	Hadlock, Hansmann, Jeanty, Yarkoni, Osaka, JSUM, Tokyo, Brenner, Doubilet, Williams
FL/AC	A ratio of Femur Length to Abdominal Circumference	Hadlock
FL/BPD	A ratio of Femur Length to Biparietal Diameter	Hohler
CI	Cephalic Index	Hadlock, Chitty
Cerebellum/AC	Transcerebellar Diameter to Abdominal Circumference Ratio	Meyer, Nicolaides
GSD	Gestational Sac, Maximum Length	Hellman, Rempen, Tokyo
LV-A	Anterior Horn Lateral ventricular width	Nicolaides
LV-P	Posterior Horn Lateral ventricular width	Nicolaides
APTD x TTD	Anteroposterior Trunk Diameter multiplied by Transverse Trunk Diameter	JSUM, Tokyo
MCA PI	Middle Cerebral Artery PI	JSUM
MCA RI	Middle Cerebral Artery RI	JSUM
Umb A PI	Umbilical Artery PI	JSUM
Umb A RI	Umbilical Artery RI	JSUM
HC/AC	Ratio of Head Circumference to Abdominal Circumference	Nicolaides, Cambell
Kidney AP	Kidney Anteroposterior	Bertagnoli, CFEF, Hansmann
Kidney Length	Kidney Length	Bertagnoli, CFEF, Hansmann
Kidney Width	Kidney Width	CFEF
MAD	Mean Abdominal Diameter	Kurmanavicius
TAD	Transverse Abdominal Diameter	CFEF
Fibula	Fibula	Chitty
Radius	Radius	Chitty
Uterine A	Uterine Artery	Merz, Schaffer

Appendix D Cardiac Reference

Cardiac Measurement and Calculation Labels	2
2D-Mode Measurement Labels	2
2D-Mode Calculation Labels.....	6
Left Ventricular Function Assessment Measurement and Calculation Labels.....	14
M-Mode Measurement Labels	16
M-Mode Calculation Labels	18
Doppler Measurement Labels.....	21
Doppler Calculation Labels.....	25
Doppler Tissue Imaging Measurement and Calculation Labels.....	30
Sections in the Patient Report.....	31

Cardiac Measurement and Calculation Labels

2D-Mode Measurement Labels

Measurement Label	Description	Measurement Method
Mitral Valve	Mitral Valve Function	---
EPSS	E Point to Septal Separation. The distance between the anterior leaflet and the ventricular septum at early diastole.	Distance
MV Area	Mitral Valve Area	Trace
MV Diam	Mitral Valve Diameter	Distance
HR	Heart Rate	HR
HR (Edit)	Heart Rate	(user-entered value)
AO/LA	Aorta/Left Atrium	---
RV Diam	Right Ventricular Diameter	Distance
AOd	Aorta at end-diastole	Distance
AOs	Aorta at end-systole	Distance
ACS	Aortic Cusp Separation	Distance
LAd Diam	Left Atrial Diameter at end-diastole	Distance
LAs Diam	Left Atrial Diameter at end-systole	Distance
Aortic Valve	Aortic Valve Function	---
LVOT Diam	Left Ventricular Outflow Tract Diameter	Distance
AV Area	Aortic Valve Area	Trace
HR	Heart Rate	HR
HR (Edit)	Heart Rate	(user-entered value)
Pulmonary Valve	Pulmonary Valve Function	---
PV Diam	Pulmonary Valve diameter	Distance
HR	Heart Rate	HR
HR Edit	Heart Rate	(user-entered value)
RV	Right Ventricle Function	---
RVOT Diam	Right Ventricular Outflow Tract Diameter	Distance
RV Diam	Right Ventricle diameter	Distance
RV Basal Diam	Right Ventricle Basal Diameter	Distance
RV Mid Diam	Right Ventricle Mid Diameter	Distance
RV Length	Right Ventricle Length	Distance
RV Area d A4C	Right Ventricular Area Diastole, Apical Four-chamber View	Trace
RV Area s A4C	Right Ventricular Area Systole, Apical Four-chamber View	Trace

Measurement Label	Description	Measurement Method
PISA (AR)	Proximal Isovelocity Surface Area (Aortic Regurgitation)	---
AR Radius	Aortic Regurgitation Radius	Distance
AR Aliasing Vel	Aortic Regurgitation Aliasing Velocity	Velocity
AR RF (Doppler)	Aortic Regurgitant Fraction by Doppler	---
AV Site 2 Diam	Aortic Valve Site 2 Diameter	Distance
LVOT Diam	Left Ventricular Outflow Tract Diameter	Distance
PISA (TR)	Proximal Isovelocity Surface Area (Tricuspid Regurgitation)	---
TR Radius	Tricuspid Regurgitation Radius	Distance
TR Aliasing Vel	Tricuspid Regurgitation Aliasing Velocity	Velocity
PISA (MR)	Proximal Isovelocity Surface Area (Mitral Regurgitation)	---
MR Radius	Mitral Regurgitation Radius	Distance
MR Aliasing Vel	Mitral Regurgitation Aliasing Velocity	Velocity
MR RF (Doppler)	Mitral Valve Regurgitant Fraction by Doppler	---
MV Diam	Mitral Valve Diameter	Distance
MV Site 2 Diam	Mitral Valve Site 2 Diameter	Distance
PISA (MS)	Proximal Isovelocity Surface Area (Mitral Stenosis)	---
MS Radius	Mitral Stenosis Radius	Distance
MS Aliasing Vel	Mitral Stenosis Aliasing Velocity	Velocity
MS Angle	Mitral Stenosis Angle	Angle
Valves	Valves	---
LVOT Diam	Left Ventricular Outflow Tract Diameter	Distance
MV Diam	Mitral Valve Diameter	Distance
TV Diam	Tricuspid Valve Diameter	Distance
PV Diam	Pulmonary Valve Diameter	Distance
MV Area	Mitral Valve Area	Spline
AV Area	Aortic Valve Area	Spline
Qp/Qs	Pulmonary-Systemic Flow Ratio	---
Systemic Diam	Systemic Diameter	Distance
Pulmonic Diam	Pulmonic Diameter	Distance
HR	Heart Rate	HR
HR (Edit)	Heart Rate	(user-entered value)

Measurement Label	Description	Measurement Method
(LV Dimensions)	Left Ventricular Dimensions	---
RVAWd	Right Ventricular Anterior Wall at end-diastole	Distance
RVDd	Right Ventricular Dimension at end-diastole	Distance
IVSd	Interventricular Septal Dimension at end-diastole	Distance
LVIDd	Left Ventricular Internal Dimension at end-diastole	Distance
LVPWd	Left Ventricular Posterior Wall Dimension at end-diastole	Distance
Systole	Left Ventricular Dimension at end-systole	Distance
Diastole	Left Ventricular Dimension at end-diastole	Distance
IVSs	Interventricular Septal Dimension at end-systole	Distance
LVIDs	Left Ventricular Internal Dimension at end-systole	Distance
LVPWs	Left Ventricular Posterior Wall Dimension at end-systole	Distance
(LV Volume-specific)	Left Ventricular Volume	---
LVA d sxPM	Left Ventricular Area at Papillary Muscle level at end-diastole in short axis view	Spline
LVA d sxMV	Left Ventricular Area at Mitral Valve level at end-diastole in short axis view	Spline
LVL d apical	Left Ventricular Length at end-diastole in apical view	Distance
LVA s sxPM	Left Ventricular Area at Papillary Muscle level at end-systole in short axis view	Spline
LVA s sxMV	Left Ventricular Area at Mitral Valve level at end-systole in short axis view	Spline
LVL s apical	Left Ventricular Length at end-systole in apical view	Distance
LVA d apical	Left Ventricular Long Axis Area at end-diastole in apical view	Spline
LVA s apical	Left Ventricular Long Axis Area at end-systole in apical view	Spline
(LV Mass-specific)	Left Ventricular Mass	---
A Sax Epi	Area short axis epicardial	Spline
A Sax Endo	Area short axis endocardial	Spline
a	Semi-major axis from widest minor axis radius to apex	Distance
d	Truncated semi-major axis from widest minor axis radius to mitral annulus plane	Distance
LVL	Left Ventricular Length	Distance

Measurement Label	Description	Measurement Method
(LA Volume-specific)	Left Atrial Volume	---
A1-A4C	End-Diastolic Volume for 4-chamber	Spline
A2-A2C	Two-chamber or apical long area of left atrium	Spline
L	Length	Distance
4CH	Four-chamber volume	Volume (Disk)
2CH	Two-chamber volume	Volume (Disk)
(RA Volume-specific)	Right Atrial Volume	---
A1-A4C	Four-chamber area of right atrium	Spline
A2-A2C	Two-chamber or apical long area of right atrium	Spline
L	Length	Distance
4CH	Four-chamber volume	Volume (Disk)
2CH	Two-chamber volume	Volume (Disk)
IVC	Inferior Vena Cava	Distance

2D-Mode Calculation Labels

Calculation Label	Description	Required Measurements	Units
CI	Cardiac Index is the cardiac output of the left ventricle per square meter of body surface area (BSA): $CI = CO \div BSA$ For height in centimeters and weight in kilograms: $BSA = 0.007184 \times (\text{Weight}^{0.425}) \times (\text{Height}^{0.725})$ For height in feet/inches and weight in pounds: $BSA = 0.007184 \times (\text{Weight} \times 0.454)^{0.425} \times (\text{Height} \times 2.54)^{0.725}$	HR and BSA (BSA is determined by height and weight of patient; entered in the patient data form) EDV (BP) ESV (BP)	L/min/m ²
LA CI	Left Atrium Cardiac Index is the cardiac output of the left atrium per square meter of body surface area (BSA): $LA\ CI = CO \div BSA$ For height in centimeters and weight in kilograms: $BSA = 0.007184 \times (\text{Weight}^{0.425}) \times (\text{Height}^{0.725})$ For height in feet/inches and weight in pounds: $BSA = 0.007184 \times (\text{Weight} \times 0.454)^{0.425} \times (\text{Height} \times 2.54)^{0.725}$	HR and BSA (BSA is determined by height and weight of patient; entered in the patient data form) LA EDV (BP) LA ESV (BP)	L/min/m ²
CI2	Cardiac Index is the cardiac output of the left ventricle per square meter of body surface area (BSA): $CI2 = CO2 \div BSA$ For height in centimeters and weight in kilograms: $BSA = 0.007184 \times (\text{Weight}^{0.425}) \times (\text{Height}^{0.725})$ For height in feet/inches and weight in pounds: $BSA = 0.007184 \times (\text{Weight} \times 0.454)^{0.425} \times (\text{Height} \times 2.54)^{0.725}$	HR and BSA (BSA is determined by height and weight of patient; entered in the patient data form) 2ch EDV 2ch ESV	L/min/m ²
LA CI2	Left Atrium Cardiac Index is the cardiac output of the left atrium per square meter of body surface area (BSA): $LA\ CI2 = LA\ CO2 \div BSA$ For height in centimeters and weight in kilograms: $BSA = 0.007184 \times (\text{Weight}^{0.425}) \times (\text{Height}^{0.725})$ For height in feet/inches and weight in pounds: $BSA = 0.007184 \times (\text{Weight} \times 0.454)^{0.425} \times (\text{Height} \times 2.54)^{0.725}$	HR and BSA (BSA is determined by height and weight of patient; entered in the patient data form) 2ch EDV 2ch ESV	L/min/m ²

Calculation Label	Description	Required Measurements	Units
CI4	Cardiac Index is the cardiac output of the left ventricle per square meter of body surface area (BSA): $CI4 = CO4 \div BSA$ For height in centimeters and weight in kilograms: $BSA = 0.007184 \times (\text{Weight}^{0.425}) \times (\text{Height}^{0.725})$ For height in feet/inches and weight in pounds: $BSA = 0.007184 \times (\text{Weight} \times 0.454)^{0.425} \times (\text{Height} \times 2.54)^{0.725}$	HR and BSA (BSA is determined by height and weight of patient; entered in the patient data form) 4ch EDV 4ch ESV	L/min/m ²
LA CI4	Left Atrium Cardiac Index is the cardiac output per of the left atrium per square meter of body surface area (BSA): $LA\ CI4 = LA\ CO4 \div BSA$ For height in centimeters and weight in kilograms: $BSA = 0.007184 \times (\text{Weight}^{0.425}) \times (\text{Height}^{0.725})$ For height in feet/inches and weight in pounds: $BSA = 0.007184 \times (\text{Weight} \times 0.454)^{0.425} \times (\text{Height} \times 2.54)^{0.725}$	HR and BSA (BSA is determined by height and weight of patient; entered in the patient data form) 4ch EDV 4ch ESV	L/min/m ²
CO	Cardiac Output is the effective volume of blood ejected from the left ventricle of the heart per unit of time: $CO = [(EDV - ESV) \div 1000](HR)$	EDV ESV HR	L/min
CO2	Cardiac Output is the effective volume of blood ejected from the left ventricle of the heart per unit of time: $CO2 = [(2ch\ EDV - 2ch\ ESV) \div 1000](HR)$	2ch EDV 2ch ESV HR	L/min
LA CO2	Left Atrium Cardiac Output is the effective volume of blood ejected from the left atrium of the heart per unit of time: $LA\ CO2 = [(2ch\ EDV - 2ch\ ESV) \div 1000](HR)$	2ch EDV 2ch ESV HR	L/min
CO4	Cardiac Output is the effective volume of blood ejected from the left ventricle of the heart per unit of time: $CO4 = [(4ch\ EDV - 4ch\ ESV) \div 1000](HR)$	4ch EDV 4ch ESV HR	L/min
LA CO4	Left Atrium Cardiac Output is the effective volume of blood ejected from the left atrium of the heart per unit of time: $LA\ CO4 = [(4ch\ EDV - 4ch\ ESV) \div 1000](HR)$	4ch EDV 4ch ESV HR	L/min

Calculation Label	Description	Required Measurements	Units
diffD(2ch)	Left Ventricle Diastolic differentiation for 2-chamber $\text{diffD}(2\text{ch}) = \text{abs}(\text{L}(\text{d})2 - \text{L}(\text{d})4) / \text{max}(\text{L}(\text{d})4, \text{L}(\text{d})2) * 100$	L(d)4 - (mm) L(d)2 - (mm)	%
diffS(2ch)	Left Ventricle Systolic differentiation for 2-chamber $\text{diffS}(2\text{ch}) = \text{abs}(\text{L}(\text{s})2 - \text{L}(\text{s})4) / \text{max}(\text{L}(\text{s})4, \text{L}(\text{s})2) * 100$	L(s)4 - (mm) L(s)2 - (mm)	%
diffD(4ch)	Left Ventricle Diastolic differentiation for 4-chamber $\text{diffD}(4\text{ch}) = \text{abs}(\text{L}(\text{d})2 - \text{L}(\text{d})4) / \text{max}(\text{L}(\text{d})4, \text{L}(\text{d})2) * 100$	L(d)4 - (mm) L(d)2 - (mm)	%
diffS(4ch)	Left Ventricle Systolic differentiation for 4-chamber $\text{diffS}(4\text{ch}) = \text{abs}(\text{L}(\text{s})2 - \text{L}(\text{s})4) / \text{max}(\text{L}(\text{s})4, \text{L}(\text{s})2) * 100$	L(s)4 - (mm) L(s)2 - (mm)	%
LA diffD(2ch)	Left Atrium Diastolic differentiation for 2-chamber $\text{LA diffD}(2\text{ch}) = \text{abs}(\text{L}(\text{d})2 - \text{L}(\text{d})4) / \text{max}(\text{L}(\text{d})4, \text{L}(\text{d})2) * 100$	L(d)4 - (mm) L(d)2 - (mm)	%
LA diffS(2ch)	Left Atrium Systolic differentiation for 2-chamber $\text{LA diffS}(2\text{ch}) = \text{abs}(\text{L}(\text{s})2 - \text{L}(\text{s})4) / \text{max}(\text{L}(\text{s})4, \text{L}(\text{s})2) * 100$	L(s)4 - (mm) L(s)2 - (mm)	%
LA diffD(4ch)	Left Atrium Diastolic differentiation for 4-chamber $\text{LA diffD}(4\text{ch}) = \text{abs}(\text{L}(\text{d})2 - \text{L}(\text{d})4) / \text{max}(\text{L}(\text{d})4, \text{L}(\text{d})2) * 100$	L(d)4 - (mm) L(d)2 - (mm)	%
LA diffS(4ch)	Left Atrium Systolic differentiation for 4-chamber $\text{LA diffS}(4\text{ch}) = \text{abs}(\text{L}(\text{s})2 - \text{L}(\text{s})4) / \text{max}(\text{L}(\text{s})4, \text{L}(\text{s})2) * 100$	L(s)4 - (mm) L(s)2 - (mm)	%
EDV (BP)	Left Ventricular End-Diastolic Volume	Varies with the selected Volume method	mL
ESV (BP)	Left Ventricular End-Systolic Volume	Varies with the selected Volume method	mL
EF (BP)	Ejection Fraction is the ratio of the stroke volume to the end-diastolic volume: $\text{EF (BP)} = 100[(\text{EDV (BP)} - \text{ESV (BP)}) \div \text{EDV}]$	EDV (BP) ESV (BP)	%
LA EDV (BP)	Left Atrium End-Diastolic Volume	Varies with the selected Volume method	mL
LA ESV (BP)	Left Atrium End-Systolic Volume	Varies with the selected Volume method	mL
LA EF (BP)	Left Atrium Ejection Fraction is the ratio of the stroke volume to the end-diastolic volume: $\text{LA EF (BP)} = 100[(\text{LA EDV (BP)} - \text{ESV}) \div \text{LA EDV (BP)}]$	LA EDV (BP) LA ESV (BP)	%

Calculation Label	Description	Required Measurements	Units
EF2	Ejection Fraction is the ratio of the stroke volume to the end-diastolic volume: $EF2 = 100[(2ch\ EDV - 2ch\ ESV) \div 2ch\ EDV]$	2ch EDV 2ch ESV	%
LA EF2	Left Atrium Ejection Fraction is the ratio of the stroke volume for 2-chamber to the end-diastolic volume: $LA\ EF2 = 100[(2ch\ ESV - 2ch\ EDV) \div 2ch\ ESV]$	LA 2ch EDV LA 2ch ESV	%
EF4	Ejection Fraction is the ratio of the stroke volume to the end-diastolic volume: $EF4 = 100[(4ch\ EDV - ESV) \div EDV]$	4ch EDV 4ch ESV	%
LA EF4	Left Atrium Ejection Fraction is the ratio of the stroke volume for 4-chamber to the end-diastolic volume: $LA\ EF4 = 100[(LA\ 4ch\ ESV - LA\ 4ch\ EDV) \div LA\ 4ch\ ESV]$	LA 4ch EDV LA 4ch ESV	%
Flow Rate	Flow Rate $Flow\ Rate = 2 * 3.1415926535 * ((Radius)^2) * Aliasing\ Vel$	Radius Aliasing Vel	mL/s
FS	Fractional Shortening is the percentage of shortening of the Left Ventricular dimension: $FS = 100[(LVIDd - LVIDs) \div LVIDd]$	LVIDd LVIDs	%
SI	Stroke Index is a measurement of left ventricle stroke volume normalized to BSA: $SI = SV \div BSA$	SV BSA	mL/m ²
LA SI	Left Atrium Stroke Index is a measurement of left atrium stroke volume normalized to BSA: $LA\ SI = LA\ SV \div BSA$	LA SV BSA	mL/m ²
SI2	Stroke Index is a measurement of left ventricle stroke volume for 2-chamber normalized to BSA: $SI2 = SV2 \div BSA$	SV2 BSA	mL/m ²
LA SI2	Left Atrium Stroke Index is a measurement of left atrium stroke volume for 2-chamber normalized to BSA: $LA\ SI2 = LA\ SV2 \div BSA$	LA SV2 BSA	mL/m ²
SI4	Stroke Index is a measurement of left ventricle stroke volume for 4-chamber normalized to BSA: $SI4 = SV4 \div BSA$	SV4 BSA	mL/m ²
LA SI4	Left Atrium Stroke Index is a measurement of left atrium stroke volume for 4-chamber normalized to BSA: $LA\ SI4 = LA\ SV4 \div BSA$	LA SV4 BSA	mL/m ²

Calculation Label	Description	Required Measurements	Units
SV	Stroke Volume is the volume of blood ejected from a ventricle during one cardiac cycle or phase of ventricular systole: $SV = EDV - ESV$	EDV (BP) ESV (BP)	mL
LA SV	Left Atrium Stroke Volume is the volume of blood ejected from an atrium during one cardiac cycle or phase of atrial systole: $LA SV = EDV (BP) - ESV (BP)$	LA EDV (BP) LA ESV (BP)	mL
SV2	Stroke Volume is the volume of blood ejected from a ventricle during one cardiac cycle or phase of ventricular systole: $SV = 2ch EDV - 2ch ESV$	2ch EDV 2ch ESV	mL
LA SV2	Left Atrium Stroke Volume is the volume of blood for 2-chamber ejected from an atrium during one cardiac cycle or phase of atrial systole: $LA SV2 = LA 2ch EDV - LA 2ch ESV$	LA 2ch EDV LA 2ch ESV	mL
SV4	Stroke Volume is the volume of blood ejected from a ventricle during one cardiac cycle or phase of ventricular systole: $SV = 4ch EDV - 4ch ESV$	4ch EDV 4ch ESV	mL
LA SV4	Left Atrium Stroke Volume is the volume of blood for 4-chamber ejected from an atrium during one cardiac cycle or phase of atrial systole: $LA SV4 = LA 4ch EDV - LA 4ch ESV$	LA 4ch EDV LA 4ch ESV	mL
t	Thickness of the myocardium $t(mm) = \sqrt{(A Sax Epi(cm^2) * 100 / 3.14) - (A Sax Endo(cm^2) * 100 / 3.14)}$	A Sax Epi A Sax Endo	mm
b	Short axis radius $b(mm) = \sqrt{(A Sax Endo(cm^2) * 100 / 3.14)}$	A Sax Endo	mm
TR Flow Rate	Tricuspid regurgitation flow rate $TR Flow Rate = 2 * 3.1415926535 * ((TR Radius)^2) * TR Aliasing Vel$	TR Radius TR Aliasing Vel	mL/s
TS Flow Rate	Tricuspid stenosis flow rate $TS Flow Rate = 2 * 3.1415926535 * ((TS Radius)^2) * TS Aliasing Vel$	TS Radius TS Aliasing Vel	mL/s

Calculation Label	Description	Required Measurements	Units
LV Mass	Left Ventricular Mass $\text{LV Mass(g)} = 1.05 * \pi * (((b + t)^2) * (2 / 3 * (a + t) + d - ((d^3) / (3 * (a + t^2)))) - ((b^2) * (2 / 3 * a + d - ((d^3) / (3 * (a^2))))))$	a d b t	g
LV Mass-I	LV Mass Index describes the proportion of left ventricular mass to the body surface area: $\text{LV Mass-I(g / m}^2\text{)} = \text{LV Mass(g)} / \text{BSA(m}^2\text{)}$	LV Mass BSA	g/m ²
LV Mass T-E	LV Mass T-E estimates left ventricular mass using a Truncated Ellipse: $\text{LV Mass(g)} = 1.05 * \pi * ((b(\text{cm}) + t(\text{cm}))^2 * (2 / 3 * (a(\text{cm}) + t(\text{cm})) + d(\text{cm}) - ((d(\text{cm}))^3 / (3 * (a(\text{cm}) + t(\text{cm}))^2))) - (b(\text{cm}))^2 * (2 / 3 * a(\text{cm}) + d(\text{cm}) - ((d(\text{cm}))^3 / (3 * (a(\text{cm}))^2))))$	a d b t	g
LV Mass A-L	LV Mass A-L estimates left ventricular volume using Area and Length: $\text{LV Mass(g)} = 1.05 * (((5 / 6) * \text{A Sax Epi(cm}^2\text{)} * (\text{LVL(cm)} + t(\text{cm}))) - (5 / 6) * \text{A Sax Endo(cm}^2\text{)} * \text{LVL(cm)}))$	A Sax Epi A Sax Endo LVL	g
LA Vol-I	LA Vol-I describes the proportion of left atrial volume to the body surface area: $\text{LA Vol-I(mL / m}^2\text{)} = \text{LA Vol(mL)} / \text{BSA(m}^2\text{)}$	LA Vol (A-L) or LA Vol (Simp) BSA	mL/m ²
LA Vol (A-L)	LA Vol (A-L) estimates left atrial volume using Area and Length: $\text{LA Vol (A-L)(mL)} = (8 / 3 / \pi) * \text{A1-A4C(cm}^2\text{)} * \text{A2-A2C(cm}^2\text{)} / \text{L(mm)} * 10$	A1-A4C A2-A2C L	mL
LA Vol (Simp)	LA Vol (Simp) uses Simpson's formula to estimate left atrial volume using the 4CH and 2CH measurements: $\text{LA Vol (Simp)(mL)} = (4\text{CH(mL)} + 2\text{CH(mL)}) / 2$	4CH 2CH	mL
LA s/Ao d	Left Atrium Systole to Aorta Diastole ratio $\text{LA s / Ao d} = \text{LA s Diam} / \text{Ao d}$	LA s Diam Ao d	---
AR Flow Rate	Aortic Regurgitation Flow Rate $\text{AR Flow Rate(mL/s)} = 2 * \pi * ((\text{AR Radius(cm)})^2) * \text{AR Aliasing Vel(cm/s)}$	AR Radius AR Aliasing Vel	mL/s
MR Flow Rate	Mitral Regurgitation flow rate $\text{MR Flow Rate} = 2 * 3.1415926535 * ((\text{MR Radius})^2) * \text{MR Aliasing Vel}$	MR Radius MR Aliasing Vel	mL/s
MS Flow Rate	Mitral Stenosis flow rate $\text{MS Flow Rate} = 2 * 3.1415926535 * ((\text{MS Radius})^2) * \text{MS Aliasing Vel}$	MS Radius MS Aliasing Vel	mL/s

Calculation Label	Description	Required Measurements	Units
RA Vol-I	RA Vol-I describes the proportion of right atrial volume to the body surface area: $\text{RA Vol-I}(\text{mL} / \text{m}^2) = \text{RA Vol}(\text{mL}) / \text{BSA}(\text{m}^2)$	RA Vol (A-L) or RA Vol (Simp) BSA	mL/m ²
RA Vol (A-L)	RA Vol (A-L) estimates right atrial volume: using Area and Length $\text{RA Vol (A-L)}(\text{mL}) = (8 / 3 / \pi) * \text{A1-A4C}(\text{cm}^2) * \text{A2-A2C}(\text{cm}^2) / \text{L}(\text{mm}) * 10$	A1-A4C A2-A2C L	mL
RA Vol (Simp)	RA Vol (Simp) uses Simpson's formula to estimate right atrial volume using the 4CH and 2CH measurements: $\text{RA Vol (Simp)}(\text{mL}) = (4\text{CH}(\text{mL}) + 2\text{CH}(\text{mL})) / 2$	4CH 2CH	mL
RA Vol-I (Simp)	Right Atrium Volume Index (Simpson) $\text{RA Vol-1 (Simp)} = \text{RA Vol (Simp)} / \text{BSA}$	RA Vol (Simp) BSA	mL/m ²
RA Vol-I (A-L)	Right Atrium Volume Index using area and length $\text{RA Vol-1 (A-L)} = \text{RA Vol (A-L)} / \text{BSA}$	RA Vol (A-L) BSA	mL/m ²
RWT	Relative Wall Thickness $\text{RWT} = 2 * \text{LVPWd} / \text{LVIDd}$	LVIDd LVPWd	mm
AOd/LAs	Aorta Diastole to Left Atrium Systole ratio	LAs AOd	---
LAs/AOd	Left Atrium Systole to Aorta Diastole ratio	LAs AOd	---
mVcfc*	Mean Velocity of Circumferential Fiber Shortening $\frac{(\text{LVIDd}(\text{mm}) - \text{LVIDs}(\text{mm}))}{\text{LVIDd}(\text{mm}) \times \text{LVETc}}$ $\text{LVETc} = \frac{(\text{LVET}(\text{ms})/1000)}{\sqrt{\text{RR}}}$ $\text{RR}(\text{sec}) = \frac{60}{\text{HR}}$	LVIDd LVIDs LVET HR	circ/s
ESWS*	End-Systolic Wall Stress $0.334 * \text{sysBP}(\text{mmHg}) * \text{LVIDs}(\text{mm}) / (\text{LVPWs}(\text{mm}) * (1 + \text{LVPWs}(\text{mm}) / \text{LVIDs}(\text{mm})))$	LVIDs LVPWs sysBP	g/cm ²

Calculation Label	Description	Required Measurements	Units
TV SV	Tricuspid Valve Stroke Volume $0.785 * ((\text{TV Diam(mm)})^2) * \text{TV VTI(cm)} / 100$	TV Diam TV VTI	mL
AR RF (PISA)	Aortic Regurgitant Fraction by Proximal Isovelocity Surface Area $(\text{AR Flow Vol(mL)} / \text{LVOT SV(mL)}) * 100$	AR Flow Vol LVOT SV	%
TR RF (PISA)	Tricuspid Valve Regurgitant Fraction by Proximal Isovelocity Surface Area $(\text{TR Flow Vol(mL)} / \text{SV(mL)}) * 100$	TR Flow Vol SV	%
MR RF (PISA)	Mitral Valve Regurgitant Fraction by Proximal Isovelocity Surface Area $(\text{MR Flow Vol(mL)} / \text{MV SV(mL)}) * 100$	MR Flow Vol MV SV	%
AR RF (Doppler)	Aortic Regurgitant Fraction by Doppler $((\text{LVOT SV(mL)} - \text{AV Site 2 SV(mL)}) / \text{LVOT SV(mL)}) * 100$	LVOT SV AV Site 2 SV	%
MR RF (Doppler)	Mitral Valve Regurgitant Fraction by Doppler $((\text{MV SV(mL)} - \text{MV Site 2 SV(mL)}) / \text{MV SV(mL)}) * 100$	MV SV MV Site 2 SV	%
AV Site 2 SV	Aortic Valve Site 2 Stroke Volume $0.785 * ((\text{AV Site 2 Diam(cm)})^2) * \text{AV Site 2 VTI(cm)}$	AV Site 2 Diam AV Site 2 VTI	mL
MV Site 2 SV	Mitral Valve Site 2 Stroke Volume $0.785 * ((\text{MV Site 2 Diam(cm)})^2) * \text{MV Site 2 VTI(cm)}$	MV Site 2 Diam MV Site 2 VTI	mL
RV FAC (%)	Right Ventricular Fractional Area Change (%) $((\text{RV Area d A4C(cm}^2\text{)} - \text{RV Area s A4C(cm}^2\text{)}) / \text{RV Area d A4C(cm}^2\text{)}) * 100$	RV Area d A4C RV Area s A4C	%
MR Volume (Doppler)	Mitral Valve Volume by Doppler $\text{MV SV(mL)} - \text{MV Site 2 SV(mL)}$	MV SV MV Site 2 SV	mL
AR Volume (Doppler)	Aortic Volume by Doppler $\text{LVOT SV(mL)} - \text{AV Site 2 SV(mL)}$	LVOT SV AV Site 2 SV	mL

* Applies only to Ped Echo exam types.

Left Ventricular Function Assessment Measurement and Calculation Labels

Volume Method	Menu Label	Description	Calculated Results	
Simpson Single Plane	Simpson SP	Estimates Volume using the Disk measurement at end-diastole and at end-systole for the apical two-chamber or four-chamber view.	EDV	mL
			ESV	mL
			SV	mL
			CO	L/min
			EF	%
			SI	mL/m ²
			CI	mL/min/m ²
Simpson Bi-Plane	Simpson BP	Estimates Volume using the Disk measurement at end-diastole and at end-systole for the apical four-chamber view: $\text{EDV (mL)} = \pi * \text{LVLd (mm)} / 20 * \sum (r_{2i}(\text{mm}) * r_{4i}(\text{mm})) / 1000$ $\text{ESV (mL)} = \pi * \text{LVLS (mm)} / 20 * \sum (r_{2i}(\text{mm}) * r_{4i}(\text{mm})) / 1000$	EDV	mL
			ESV	mL
			SV	mL
			CO	L/min
			EF	%
			SI	mL/m ²
			CI	mL/min/m ²
Modified Simpson Formula	Mod. Simpson	Estimates Volume using the trace measurement in two different views at end-diastole and at end-systole for sax MV and for sax PM. Also requires a distance measurement of the long axis dimension in an apical four-chamber view at end-diastole and at end-systole: $\text{EDV (mL)} = \text{LVLd apical (mm)} / 9 * (4 * \text{LVAd sax MV (cm}^2\text{)} + 2 * \text{LVAd sax PM (cm}^2\text{)} + \text{SQRT (LVAd sax MV (cm}^2\text{)} * \text{LVAd sax PM (cm}^2\text{)})} / 10)$ $\text{ESV (mL)} = \text{LVLS apical (mm)} / 9 * (4 * \text{LVAs sax MV (cm}^2\text{)} + 2 * \text{LVAs sax PM (cm}^2\text{)} + \text{SQRT (LVAs sax MV (cm}^2\text{)} * \text{LVAs sax PM (cm}^2\text{)})} / 10)$	EDV	mL
			ESV	mL
			SV	mL
			CO	L/min
			EF	%
			SI	mL/m ²
			CI	mL/min/m ²
Single Plane Formula	Single Plane	Estimates Volume using the trace measurement in a single plane at end-diastole and at end-systole for apical two- or four-chamber view. Also requires a distance measurement of the long axis dimension in an apical view at end-diastole end-systole: $\text{EDV} = (8 \div (3\pi)) * (\text{LVAd apical (cm}^2\text{)})^2 / \text{LVLd apical (mm)} * 10$ $\text{ESV} = (8 \div (3\pi)) * (\text{LVAs apical (cm}^2\text{)})^2 / \text{LVLS apical (mm)} * 10$	EDV	mL
			ESV	mL
			SV	mL
			CO	L/min
			EF	%
			SI	mL/m ²
			CI	mL/min/m ²

Volume Method	Menu Label	Description	Calculated Results	
Bi-plane Formula	Bi-Plane	Estimates Volume using the trace measurement in two planes at end-diastole and at end-systole for apical and for sax MV. Also requires a distance measurement of the LVIDd and the LVIDs: $EDV = (8 \div (3\pi)) * LVAd\ apical(cm^2) * LVAd\ sax\ MV(cm^2) / LVIDd(mm) * 10$ $ESV = (8 \div (3\pi)) * LVAs\ apical(cm^2) * LVAs\ sax\ MV(cm^2) / LVIDs(mm) * 10$	EDV	mL
			ESV	mL
			SV	mL
			CO	L/min
			EF	%
			SI	mL/m ²
			CI	mL/min/m ²
Bullet Formula	Bullet	Estimates Volume using the trace measurement in one view at end-diastole and at end-systole for sax MV. Also requires a distance measurement of the long axis dimension in an apical view at end-diastole end-systole. This formula assumes that the left ventricle is shaped like a bullet; the base being cylindrical and the apex cone-like: $EDV = (5 \div 6) * LVLd\ apical(mm) * LVAd\ sax\ MV(cm^2) / 10$ $ESV = (5 \div 6) * LVLs\ apical(mm) * LVAs\ sax\ MV(cm^2) / 10$	EDV	mL
			ESV	mL
			SV	mL
			CO	L/min
			EF	%
			SI	mL/m ²
			CI	mL/min/m ²
Cubed Formula	Cubed	Estimates Volume from the Left Ventricular Internal Diameter (LVID) using the distance measurement: $EDV = LVIDd(mm)^3 / 1000$ $ESV = LVIDs(mm)^3 / 1000$	EDV	mL
			ESV	mL
			SV	mL
			CO	L/min
			EF	%
			SI	mL/m ²
			CI	mL/min/m ²
Teichholz Formula	Teichholz	Estimates Volume from the Left Ventricular Internal Diameter (LVID) using the distance measurement: $EDV(mL) = 7 * (LVIDd(mm) * LVIDd(mm) * LVIDd(mm) / 1000) / (2.4 + LVIDd(mm) / 10)$ $ESV(mL) = (7 * LVIDs(mm) * LVIDs(mm) * LVIDs(mm) / 1000) / (2.4 + LVIDs(mm) / 10)$	EDV	mL
			ESV	mL
			SV	mL
			CO	L/min
			EF	%
			SI	mL/m ²
			CI	mL/min/m ²
Gibson Formula	Gibson	Estimates Volume from the Left Ventricular Internal Diameter (LVID) using the distance measurement: $EDV(mL) = \pi / 6 * (0.98 * LVIDd(mm) / 10 + 5.90) * LVIDd(mm) / 10 * LVIDd(mm) / 10$ $ESV(mL) = \pi / 6 * (1.14 * LVIDs(mm) / 10 + 4.18) * LVIDs(mm) / 10 * LVIDs(mm) / 10$	EDV	mL
			ESV	mL
			SV	mL
			CO	L/min
			EF	%
			SI	mL/m ²
			CI	mL/min/m ²

M-Mode Measurement Labels

Measurement Label	Description	Measurement Method
AO/LA	Aortic Valve/Left Atrium function	---
AV/LA(M)	Aortic Valve/Left Atrium function	Distance
RV Diam	Right Ventricular Diameter	Distance
AOd	Aorta at end-diastole	Distance
AOs	Aorta at end-systole	Distance
ACS	Aortic Cusp Separation	Distance
LAd Diam	Left Atrial Diameter at end-diastole	Distance
LAs Diam	Left Atrial Diameter at end-systole	Distance
LV ET	Left Ventricular Ejection Time	Time
LV PEP	Left Ventricular Pre-Ejection Period	Time
Mitral Valve	Mitral Valve function	---
CE amp	Amplitude of the E wave	Distance
CA amp	Amplitude of the A wave	Distance
DE Slope	Anterior excursion of the mitral leaflet at early diastole	Slope
DE excursion	Amplitude of the early diastole wave	Distance
EPSS	E Point to Septal Separation. The distance between the anterior leaflet and the ventricular septum at early diastole.	Distance
EF Slope	Deceleration slope of the mitral leaflet's closure at early diastole	Slope
MAPSE	Mitral Annular Plane Systolic Excursion	Distance
RV	Right Ventricular Dimensions	---
RV Diam	Right Ventricular Diameter	Distance
TAPSE	Tricuspid Annular Plane Systolic Excursion	Distance

Measurement Label	Description	Measurement Method
(LV Dimensions)	Left Ventricle Dimensions	---
Diastole	Left Ventricular Dimension at end-diastole	3Distance
Systole	Left Ventricular Dimension at end-systole	3Distance
RVDd	Right Ventricular Dimension at end-diastole	Distance
IVSd	Interventricular Septal dimension at end-diastole	Distance
LVIDd	Left Ventricular Internal dimension at end-diastole	Distance
LVPWd	Left Ventricular Posterior wall dimension at end-diastole	Distance
IVSs	Interventricular Septal dimension at end-systole	Distance
LVIDs	Left Ventricular Internal Dimension at end-systole	Distance
LVPWs	Left Ventricular Posterior Wall dimension at end-systole	Distance
(LV Function)	Left Ventricular Function	---
LV ET	Left Ventricular Ejection time	Time
LV ET (Edit)	Left Ventricular Ejection time	(user-entered value)
HR	Heart Rate	HR
HR (Edit)	Heart Rate	(user-entered value)

M-Mode Calculation Labels

Calculation Label	Description	Required Measurements	Units
CI	Cardiac Index is the cardiac output of the left ventricle per square meter of body surface area (BSA): $CI = CO \div BSA$ For height in centimeters and weight in kilograms: $BSA = 0.007184 \times (\text{Weight}^{0.425}) \times (\text{Height}^{0.725})$ For height in feet/inches and weight in pounds: $BSA = 0.007184 \times (\text{Weight} \times 0.454)^{0.425} \times (\text{Height} \times 2.54)^{0.725}$	BSA entry in the patient data form (determined by height and weight of patient)	L/min/m ²
LA CI	Left Atrium Cardiac Index is the cardiac output of the left atrium per square meter of body surface area (BSA): $LA\ CI = LA\ CO \div BSA$ For height in centimeters and weight in kilograms: $BSA = 0.007184 \times (\text{Weight}^{0.425}) \times (\text{Height}^{0.725})$ For height in feet/inches and weight in pounds: $BSA = 0.007184 \times (\text{Weight} \times 0.454)^{0.425} \times (\text{Height} \times 2.54)^{0.725}$	BSA entry in the patient data form (determined by height and weight of patient)	L/min/m ²
CO	Cardiac Output is the effective volume of blood ejected by either ventricle of the heart per unit of time: $CO = [(EDV - ESV) \div 1000](HR)$	EDV ESV HR	L/min
LA CO	Left Atrium Cardiac Output is the effective volume of blood ejected by either atrium of the heart per unit of time: $LA\ CO = [(LA\ ESV\ (BP) - LA\ EDV\ (BP)) \div 1000](HR)$	LA EDV (BP) LA ESV (BP) HR	L/min
EDV and ESV	Left Ventricular End-Systolic Volume and Left Ventricular End-Diastolic Volume estimates volume from the left ventricular internal diameter (LVID) using the distance measurement. Cubed Volume Method: $EDV = LVIDd^3$ $ESV = LVIDs^3$ Teichholz Volume Method: $EDV = 7(LVIDd^3 / 1000) / (2.4 + LVIDd/10)$ $ESV = 7(LVIDs^3 / 1000) / (2.4 + LVIDs/10)$ Gibson Volume Method: $EDV = \pi / 6 * (0.98 * LVIDd/10 + 5.90) * LVIDd/10 * LVIDd/10$ $ESV = \pi / 6 * (1.14 * LVIDs/10 + 4.18) * LVIDs/10 * LVIDs/10$	LVIDd LVIDs	mL
EF	Ejection Fraction is the ratio of the stroke volume to the end-diastolic volume: $EF = 100[(EDV - ESV) \div EDV]$	EDV ESV	%

Calculation Label	Description	Required Measurements	Units
AOd/LAs	Aorta Diastole to Left Atrium Systole ratio	LAs AOd	---
LAs/AOd	Left Atrium Systole to Aorta Diastole ratio	LAs AOd	---
LV STI	Left Ventricular Systolic Time Intervals $LV\ STI = LV\ PEP \div LV\ ET$	LV PEP LV ET	---
RWT	Relative Wall Thickness	LVID LVPWd	mm
HR	Heart Rate: $HR = 60 \div (R-R\ interval)$	HR a cardiac cycle	bpm
SI	Stroke Index is a measurement of left ventricle stroke volume normalized to BSA: $SI = SV \div BSA$	SV BSA	mL/m ²
LA SI	Left Atrium Stroke Index is a measurement of left atrium stroke volume normalized to BSA: $LA\ SI = LA\ SV \div BSA$	LA SV BSA	mL/m ²
SV	Stroke Volume is the volume of blood ejected from a ventricle during one cardiac cycle: $SV = EDV - ESV$	EDV ESV	mL
LA SV	Left Atrium Stroke Volume is the volume of blood ejected from an atrium during one cardiac cycle: $LA\ SV = LA\ EDV\ (BP) - LA\ ESV\ (BP)$	LA EDV (BP) LA ESV (BP)	mL
CA/CE	Mitral Valve ratio of the CA amplitude to the CE amplitude $CA\ amp / CE\ amp$	CA amp CE amp	---

Calculation Label	Description	Required Measurements	Units
LV Mass	Left Ventricular Mass $\text{LV Mass(g)} = 1.04 * ((\text{IVSd(mm)} / 10 + \text{LVIDd(mm)} / 10 + \text{LVPWd(mm)} / 10)^3 - (\text{LVIDd(mm)} / 10)^3)$	IVSd LVIDd LVPWd	g
LV Mass-c	LV Mass-corrected $\text{LV Mass-c(g)} = 0.8 * \text{LV Mass(g)} + 0.6$	LV Mass	g
LV Mass-I	LV Mass-Index describes the proportion of left ventricular mass to the body surface area $\text{LV Mass-I(g/m}^2\text{)} = \text{LV Mass-c(g)} / \text{BSA(m}^2\text{)}$	LV Mass-c BSA	g/m ²
mVcf	M-mode Velocity of Circumferential Fiber Shortening $((\text{LVIDd(mm)} - \text{LVIDs(mm)}) / (\text{LVIDd(mm)} * \text{LVET(ms)} / 1000))$	LVIDd LVIDs LV ET	circ/s
mVcfc*	Mean Velocity of Circumferential Fiber Shortening $\frac{(\text{LVIDd(mm)} - \text{LVIDs(mm)})}{\text{LVIDd(mm)} \times \text{LVETc}}$ $\text{LVETc} = \frac{(\text{LVET(ms)} / 1000)}{\sqrt{\text{RR}}}$ $\text{RR(sec)} = \frac{60}{\text{HR}}$	LVIDd LVIDs LV ET HR	circ/s
ESWS*	End-Systolic Wall Stress $0.334 * \text{sysBP(mmHg)} * \text{LVIDs(mm)} / (\text{LVPWs(mm)} * (1 + \text{LVPWs(mm)} / \text{LVIDs(mm)}))$	LVIDs LVPWs sysBP	g/cm ²

*Applies only to Ped Echo exam.

Doppler Measurement Labels

Measurement Label	Description	Measurement Method
Aortic Valve	Aortic Valve Function	---
Ao Arch VTI	Aortic Arch Velocity-Time Integral	VTI
Ao Arch Vmax	Aortic Arch Maximum Velocity	Velocity
AV VTI	Aortic Valve Velocity-Time Integral	VTI
LVOT VTI	Left Ventricular Outflow Tract Velocity-Time Integral	VTI
IVRT	Isovolumic Relaxation Time	Time
HR	Heart Rate	HR
HR (Edit)	Heart Rate	(user-entered value)
AV VTI	Aortic Valve Velocity-Time Integral	VTI
AV Vmax	Aortic Valve Maximum Velocity	Velocity
AV Accel Time	Aortic Valve Acceleration Time	Acceleration
LVOT Vmax	Left Ventricular Outflow Tract Maximum Velocity	Velocity
VSD Vmax	Ventricular Septal Defect Maximum Velocity	Velocity
LV ET	Left Ventricular Ejection Time	Time
LV PEP	Left Ventricular Pre-Ejection Period	Time
Mitral Valve	Mitral Valve Function	---
E Dur	E-wave Duration	Time
A Dur	A-wave Duration	Time
IVRT	Isovolumic Relaxation Time	Time
MV Decel Time	Mitral Valve Deceleration Time	Acceleration
MV E/A Vel	Mitral Valve E point/A point Velocity	Velocity
MV Peak A vel	Mitral Valve A point Velocity	Velocity
MV Peak E vel	Mitral Valve E point Velocity	Velocity
MV VTI	Mitral Valve Velocity-Time Integral	VTI
MV Vmax	Mitral Valve Maximum Velocity	Velocity
HR	Heart Rate	HR
HR (Edit)	Heart Rate	(user-entered value)
IVCT	Isovolumic Contraction Time	Time
Tricuspid Valve	Tricuspid Valve Function	---
TV Vmax	Tricuspid Valve Velocity maximum	Velocity
RAP select	Right Atrial Pressure selection	ComboType
TR VTI	Tricuspid Regurgitant Velocity-Time Integral	VTI
TV E pt	Tricuspid Valve E point	Velocity
TV A pt	Tricuspid Valve A point	Velocity
TV VTI	Tricuspid Valve Velocity-Time Integral	VTI

Measurement Label	Description	Measurement Method
Pulmonic Valve	Pulmonary Valve Function	---
PR Vmax	Pulmonary Regurgitation Velocity maximum	Velocity
PV Vmax	Pulmonary Valve Velocity maximum	Velocity
PV Accel Time	Pulmonary Valve Acceleration Time	Acceleration
RV ET	Right Ventricular Ejection Time	Time
RV Accel Time	Right Ventricular Acceleration Time	Time
RV PEP	Right Ventricular Pre-ejection Period	Time
PV VTI	Pulmonary Valve, Velocity-Time Integral	VTI
RVOT VTI	Right Ventricular Outflow Tract, Velocity-Time Integral	VTI
RVOT Vmax	Right Ventricular Outflow Tract, Velocity maximum	Velocity
HR	Heart Rate	HR
Pulmonary Vein	Pulmonary Vein Function	---
PVs1 Vel	Pulmonary Vein systole1 Velocity	Velocity
PVs2 Vel	Pulmonary Vein systole2 Velocity	Velocity
PVd Vel	Pulmonary Vein diastolic Velocity	Velocity
PVa Vel	Pulmonary Vein late diastolic Velocity	Velocity
PVa dur	Pulmonary Vein late diastolic duration	Time
PVs VTI	Pulmonary Vein systolic Velocity-Time Integral	VTI
PVd VTI	Pulmonary Vein diastolic Velocity-Time Integral	VTI
PVd Decel	Pulmonary Vein diastolic Deceleration Time	Acceleration
AR	Aortic Regurgitation	---
AR Decel Time	Deceleration Time	Acceleration
AR PHT	Aortic Insufficiency Pressure Half-Time	PHT
AR Vmax	Aortic Regurgitation Velocity maximum	Velocity
TR	Tricuspid Regurgitation	---
TR Vmean	Tricuspid Regurgitation Velocity mean	VTI
TR VTI	Tricuspid Regurgitant Velocity-Time Integral	VTI
TR Vmax	Tricuspid Regurgitation Velocity maximum	Velocity
RAP select	Right Atrial Pressure selection	ComboType
RAP(edit)	Right Atrial Pressure selection	(user-entered value)
PR	Pulmonary Regurgitation	---
PR Vmean	Pulmonary Regurgitation Velocity mean	VTI
PR Vmax	Pulmonary Regurgitation Velocity maximum	Velocity
PR Ved	Pulmonary Regurgitation Velocity end diastole	Velocity
PR VTI	Pulmonary Regurgitation Velocity-Time Integral	VTI
RAP select	Right Atrial Pressure selection	ComboType

Measurement Label	Description	Measurement Method
MR	Mitral Regurgitation	---
MR Vmax	Mitral Regurgitation Velocity maximum	Velocity
dP/dt	Rate of change of pressure over time $dP / dt = 32 / dt \times 1000$	dP/dt
LV dP/dt	Rate of change of pressure over time	dP/dt
PISA (AR)	Proximal Isovelocity Surface Area (Aortic Regurgitation)	---
AR VTI	Aortic Regurgitation Velocity-Time Integral	VTI
AR RF (Doppler)	Aortic Regurgitant Fraction by Doppler	---
AV Site 2 VTI	Aortic Valve Site 2 Velocity-Time Integral	VTI
LVOT VTI	Left Ventricular Outflow Tract Velocity-Time Integral	VTI
PISA (TR)	Proximal Isovelocity Surface Area (Tricuspid Regurgitation)	---
TR VTI	Tricuspid Regurgitation Velocity-Time Integral	VTI
PISA (MR)	Proximal Isovelocity Surface Area (Mitral Regurgitation)	---
MR VTI	Mitral Regurgitation Velocity-Time Integral	VTI
MR RF (Doppler)	Mitral Valve Regurgitant Fraction by Doppler	---
MV VTI	Mitral Valve Velocity-Time Integral	VTI
MV Site 2 VTI	Mitral Valve Site 2 Velocity-Time Integral	VTI
PISA (MS)	Proximal Isovelocity Surface Area (Mitral Stenosis)	---
MS VTI	Mitral Stenosis Velocity-Time Integral	VTI
MS Vmax	Mitral Stenosis Velocity maximum	Velocity
DTI (MV)	Doppler Tissue Imaging (Mitral Valve)	---
MV Med E' Vel	Early Diastolic Myocardial Velocity, Medial in Mitral Valve	Velocity
MV Med A' Vel	Late Diastolic Myocardial Velocity, Medial in Mitral Valve	Velocity
MV Med S' Vel	Systolic Myocardial Velocity, Medial in Mitral Valve	Velocity
MV Lat E' Vel	Early Diastolic Myocardial Velocity, Lateral in Mitral Valve	Velocity
MV Lat A' Vel	Late Diastolic Myocardial Velocity, Lateral in Mitral Valve	Velocity
MV Lat S' Vel	Systolic Myocardial Velocity, Lateral in Mitral Valve	Velocity
DTI (RV)	Doppler Tissue Imaging (Right Ventricle)	---
RV E' Vel	Early diastolic motion	Velocity
RV A' Vel	Late diastolic motion	Velocity
RV S' Vel	Systolic motion	Velocity

Measurement Label	Description	Measurement Method
Tei-Index	Tei-Index	---
LV ET	Left Ventricle Ejection Time	Time
MV C-O Time	Mitral Valve Close Open duration	Time
RV ET	Right Ventricular Ejection Time	Time
TV C-O Time	Tricuspid Valve Close-Open duration	Time
Qp/Qs	Pulmonary Valve / Aortic Valve	---
Systemic VTI	Systemic Velocity-Time Integral	VTI
Pulmonic VTI	Pulmonic Velocity-Time Integral	VTI
HR	Heart Rate	HR
HR (Edit)	Heart Rate	(user-entered value)
IVC	Inferior Vena Cava	Velocity

Doppler Calculation Labels

Calculation Label	Description	Required Measurements	Units
A/E	Mitral Valve ratio of the A point to the E point $A/E = \text{MV Peak A vel (m/s)} / \text{MV Peak E vel (m/s)}$	MV Peak A vel MV Peak E vel	---
E/A	Mitral Valve ratio of the E point to the A point $E/A = \text{MV Peak E vel (m/s)} / \text{MV Peak A vel (m/s)}$	MV Peak E vel MV Peak A vel	---
E/A	Tricuspid Valve ratio of the E point to the A point $E/A = \text{TV E pt (m/s)} / \text{TV A pt (m/s)}$	TV E pt TV A pt	---
MVA (PHT)	Mitral Valve Area (Pressure Half-Time) $MVA (PHT) = 220 / \text{MV PHT}$	MV Peak E vel MV Peak A vel	cm ²
SV	Stroke Volume $SV(\text{mL}) = 0.785 * (\text{MV diam}(\text{mm}))^2 * \text{MV VTI}(\text{cm}) / 100$	MV VTI MV diam	mL
SI	Stroke Index $SI(\text{mL/m}^2) = \text{SV (mL)} / \text{BSA (m}^2\text{)}$	SV BSA	mL/m ²
CO	Cardiac Output (Doppler-derived) $CO(\text{L/min}) = \text{SV(mL)} * \text{HR(bpm)}$	SV Heart Rate	L/min
CI	Cardiac Index $CI = \text{CO(L/min)} / \text{BSA(m}^2\text{)}$	CO BSA	L/min/m ²
LV dP/dt	Rate of change of pressure over time $LV dP/dt = 32000 / dt (\text{ms})$	---	mmHg/s
AR EO Area	Aortic Regurgitation Effective Regurgitant Orifice Area $AR EO Area(\text{mm}^2) = \text{AR Flow Rate}(\text{mL/s}) / \text{AR Vmax}(\text{m/s})$	AR Flow Rate AR Vmax	mm ²
AR Flow Rate	Aortic Regurgitation Flow Rate $AR Flow Rate(\text{mL/s}) = 2 * \pi * ((\text{AR Radius}(\text{cm}))^2) * \text{AR Aliasing Vel}(\text{cm/s})$	AR Radius AR Aliasing Vel	mL/s
AR Flow Vol	Aortic Regurgitation Flow Volume $AR Flow Vol(\text{mL}) = \text{AR EO Area}(\text{cm}^2) * \text{AR VTI}(\text{cm})$	AR EO Area AR VTI	mL

Calculation Label	Description	Required Measurements	Units
MR Flow Vol	Mitral Regurgitation Flow Volume $\text{MR Flow Vol(mL)} = \text{MR EO Area(cm}^2\text{)} * \text{MR VTI(cm)}$	MR EO Area MR VTI	mL
MR Flow Rate	Mitral Regurgitation Flow Rate $\text{MR Flow Rate(mL/s)} = 2 * \pi * ((\text{MR Radius(mm)})^2) * \text{MR Aliasing Vel(m/s)}$	MR Radius MR Aliasing Vel	mL/s
MR EO Area	Mitral Regurgitation EO Area $\text{MR EO Area(mm}^2\text{)} = \text{MR Flow Rate(mL/s)} / \text{MR Vmax(m/s)}$	MR Flow Rate MR Vmax	mm ²
MS MVA	Mitral Stenosis Mitral Valve Area $\text{MS MVA(cm}^2\text{)} = 2 * \pi * ((\text{MS Radius(cm)})^2) * \text{MS Aliasing Vel(m/s)} / \text{MS Vmax(m/s)} * \text{MS Angle(}^\circ\text{)} / 180$	MS Radius MS Aliasing Vel MS Vmax MS Angle	cm ²
MS Flow Rate	Mitral Stenosis Flow Rate $\text{MS Flow Rate(mL/s)} = 2 * \pi * ((\text{MS Radius(mm)})^2) * \text{MS Aliasing Vel(m/s)}$	MS Radius MS Aliasing Vel	mL/s
MS Flow Vol	Mitral Stenosis Flow Volume $\text{MS Flow Vol(mL)} = \text{MS MVA(cm}^2\text{)} * \text{MS VTI(cm)}$	MS MVA MS VTI	mL
LV STI	Left Ventricular Systolic Time Intervals $\text{LV STI} = \text{LV PEP(ms)} / \text{LV ET(ms)}$	LV PEP LV ET	---
PAEDP	Pulmonary Artery End Diastolic Pressure $\text{PAEDP(mmHg)} = \text{RAP(mmHg)} + (\text{PR Ved(m/s)})^2 * 4$	RAP PR Ved	mmHg
PVs2/PVd	Pulmonary Vein systole2 Velocity / Pulmonary Vein diastolic Velocity $\text{PVs2/PVd} = \text{PVs2 Vel (m/s)} / \text{PVd Vel (m/s)}$	PVs2 Vel PVd Vel	---
Sys Fraction	Systolic Filling Fraction $\text{Sys Fraction} = (\text{PVs VTI(cm)} / (\text{PVs VTI(cm)} + \text{PVd VTI(cm)}))$	PVs VTI PVd VTI	---
RVSP	Right Ventricular Systolic Pressure is determined from the tricuspid regurgitation velocity converted to a pressure gradient: $\text{RVSP(mmHg)} = (\text{TR Vmax(m/s)})^2 * 4 + \text{RAP(mmHg)}$	TR Vmax RAP	mmHg
LVOT/AV VTI	Left Ventricular Outflow Tract Aortic Valve Velocity-Time Integral $\text{LVOT/AV VTI} = \text{LVOT VTI (cm)} / \text{AV VTI(cm)}$	LVOT VTI AV VTI	---
LV MPI	Left Ventricular Index of Myocardial Performance $\text{LV MPI} = (\text{MV C-O Time (ms)} - \text{LVET(ms)}) / \text{LVET(ms)}$	MV C-O Time LV ET	---
RV MPI	Right Ventricular Index of Myocardial Performance $\text{RV MPI} = (\text{TV C-O Time (ms)} - \text{RVET(ms)}) / \text{RVET(ms)}$	TV C-O Time RV ET	---

Calculation Label	Description	Required Measurements	Units
HR	Heart Rate: HR = 60 / R-R interval.	Measure across one heart cycle	bpm
E'/A' (Med)	Doppler Tissue Imaging ratio of the lateral or medial E' and A' measurements $E'/A'(\text{Med}) = \text{Med E' Vel}(\text{m/s}) / \text{Med A' Vel}(\text{m/s})$	Med E' Vel Med A' Vel	---
E'/A' (Lat)	Doppler Tissue Imaging ratio of the lateral E' and A' measurements $E'/A'(\text{Lat}) = \text{Lat E' Vel}(\text{m/s}) / \text{Lat A' Vel}(\text{m/s})$	Lat E' Vel Lat A' Vel	---
E/E' (Lat)	Doppler Tissue Imaging ratio of the lateral E and E' measurements $E/E'(\text{Lat}) = \text{MV Peak E vel}(\text{m/s}) / \text{Lat E' Vel}(\text{m/s})$	MV Peak E vel Lat E' Vel	---
E/E' (Med)	Doppler Tissue Imaging ratio of the medial E and E' measurements $E/E'(\text{Med}) = \text{MV Peak E vel}(\text{m/s}) / \text{Med E' Vel}(\text{m/s})$	MV Peak E vel Med E' Vel	---
MVA (VTI)	Mitral Valve Area (Velocity-Time Integral) $MVA(\text{VTI})(\text{cm}^2) = \pi * (\text{LVOT Diam}(\text{mm})/2/10)^2 * (\text{LVOT VTI}(\text{cm}) / \text{MV VTI}(\text{cm}))$	LVOT Diam LVOT VTI MV VTI	cm ²
AVA (VTI)	Aortic Valve Area (Velocity-Time Integral) $AVA(\text{VTI})(\text{cm}^2) = \pi * (\text{LVOT Diam}(\text{cm})/2)^2 * (\text{LVOT VTI}(\text{cm}) / \text{AV VTI}(\text{cm}))$	LVOT Diam LVOT VTI AV VTI	cm ²
AVA (Vmax)	Aortic Valve Area Maximum Velocity $AVA(\text{Vmax}) = \pi * (\text{LVOT Diam}(\text{cm})/2)^2 * (\text{LVOT Vmax}(\text{m/s}) / \text{AV Vmax}(\text{m/s}))$	LVOT Diam LVOT Vmax AV Vmax	cm ²
RV Acc T/ET	RV Acc T/ET = RV Accel Time(ms) / RV ET(ms)	RV ET RV Accel Time	---
RV STI	Right Ventricular Systolic Time Intervals RV STI = RV PEP(m/s) / RV ET(m/s)	RV PEP RV ET	---
Qp/Qs	$Qp/Qs = \text{CO}(\text{L/min})(\text{Pulmonary Valve}) / \text{CO}(\text{L/min})(\text{Aortic Valve})$	CO Pulmonary Valve CO Aortic Valve	---
SVs	SV(mL)(Aortic Valve)	SV Aortic Valve	mL
COs	CO(L/min)(Aortic Valve)	CO Aortic Valve	L/min
SVp	SV(mL)(Pulmonary Valve)	SV Pulmonary Valve	mL
COp	CO(L/min)(Pulmonary Valve)	CO Pulmonary Valve	L/min
Qp-Qs	$Qp-Qs(\text{L/min}) = \text{CO}(\text{L/min})(\text{Pulmonary Valve}) - \text{CO}(\text{L/min})(\text{Aortic Valve})$	CO Pulmonary Valve CO Aortic Valve	L/min

Calculation Label	Description	Required Measurements	Units
MR PGmax	Mitral Regurgitation Maximum Pressure Gradient $4 * ((MR Vmax)^2)$	MR Vmax	mmHg
MS PGmax	Mitral Stenosis Maximum Pressure Gradient $4 * ((MS Vmax)^2)$	MS Vmax	mmHg
TR PGmax	Tricuspid Regurgitation Maximum Pressure Gradient $4 * ((TR Vmax)^2)$	TR Vmax	mmHg
TR EO Area	Tricuspid Regurgitation EO Area TR Flow Rate / TR Vmax	TR Flow Rate TR Vmax	mm ²
MV PGmax	Mitral Valve Maximum Pressure Gradient $4 * ((MV Vmax)^2)$	MV Vmax	mmHg
TR Flow Vol	Tricuspid Regurgitation Flow Volume TR EO Area * TR VTI	TR EO Area TR VTI	mL
TR Flow Rate	Tricuspid Regurgitation Flow Rate $TR Flow Rate(mL/s) = 2 * \pi * ((TR Radius(cm))^2) * TR Aliasing Vel(cm/s)g$	TR Radius TR Aliasing Vel	mL/s
AV PGmax	Aortic Valve Maximum Pressure Gradient $4 * ((AV Vmax)^2)$	AV Vmax	mmHg
Ao Arch PGmax	Aortic Arch Maximum Pressure Gradient $4 * ((Ao Arch Vmax)^2)$	Ao Arch Vmax	mmHg
LVOT PGmax	Left Ventricular Outflow Tract Maximum Pressure Gradient $4 * ((LVOT Vmax)^2)$	LVOT Vmax	mmHg
VSD PGmax	Ventricular Septal Defect maximum Pressure Gradient $4 * ((VSD Vmax)^2)$	VSD Vmax	mmHg
AV PGmax	Aortic Valve maximum Pressure Gradient $4 * ((AV Vmax)^2)$	AV Vmax	mmHg
AR PGmax	Aortic Regurgitation maximum Pressure Gradient $4 * ((AR Vmax)^2)$	AR Vmax	mmHg
PV PGmax	Pulmonary Valve Maximum Pressure Gradient $4 * ((PV Vmax)^2)$	PV Vmax	mmHg
PR PGmax	Pulmonary Regurgitation Maximum Pressure Gradient $4 * ((PR Vmax)^2)$	PR Vmax	mmHg
TV PGmax	Tricuspid Valve Maximum Pressure Gradient $4 * ((TV Vmax)^2)$	TV Vmax	mmHg

Calculation Label	Description	Required Measurements	Units
TV SV	Tricuspid Valve Stroke Volume $0.785 * ((\text{TV Diam}(\text{mm}))^2 * \text{TV VTI}(\text{cm})) / 100$	TV Diam TV VTI	mL
AR RF (PISA)	Aortic Regurgitant Fraction by Proximal Isovelocity Surface Area $(\text{AR Flow Vol}(\text{mL}) / \text{LVOT SV}(\text{mL})) * 100$	AR Flow Vol LVOT SV	%
TR RF (PISA)	Tricuspid Valve Regurgitant Fraction by Proximal Isovelocity Surface Area $(\text{TR Flow Vol}(\text{mL}) / \text{SV}(\text{mL})) * 100$	TR Flow Vol SV	%
MR RF (PISA)	Mitral Valve Regurgitant Fraction by Proximal Isovelocity Surface Area $(\text{MR Flow Vol}(\text{mL}) / \text{MV SV}(\text{mL})) * 100$	MR Flow Vol MV SV	%
AR RF (Doppler)	Aortic Regurgitant Fraction by Doppler $((\text{LVOT SV}(\text{mL}) - \text{AV Site 2 SV}(\text{mL})) / \text{LVOT SV}(\text{mL})) * 100$	LVOT SV AV Site 2 SV	%
MR RF (Doppler)	Mitral Valve Regurgitant Fraction by Doppler $((\text{MV SV}(\text{mL}) - \text{MV Site 2 SV}(\text{mL})) / \text{MV SV}(\text{mL})) * 100$	MV SV MV Site 2 SV	%
AV Site 2 SV	Aortic Valve Site 2 Stroke Volume $0.785 * ((\text{AV Site 2 Diam}(\text{cm}))^2 * \text{AV Site 2 VTI}(\text{cm}))$	AV Site 2 Diam AV Site 2 VTI	mL
MV Site 2 SV	Mitral Valve Site 2 Stroke Volume $0.785 * ((\text{MV Site 2 Diam}(\text{cm}))^2 * \text{MV Site 2 VTI}(\text{cm}))$	MV Site 2 Diam MV Site 2 VTI	mL
RVOT PGmax	Right Ventricular Outflow Tract Maximum Pressure Gradient $4 * ((\text{RVOT Vmax}(\text{m/s}))^2)$	RVOT Vmax	mmHg
MR Volume (Doppler)	Mitral Valve Volume by Doppler $\text{MV SV}(\text{mL}) - \text{MV Site 2 SV}(\text{mL})$	MV SV MV Site 2 SV	mL
AR Volume (Doppler)	Aortic Volume by Doppler $\text{LVOT SV}(\text{mL}) - \text{AV Site 2 SV}(\text{mL})$	LVOT SV AV Site 2 SV	mL

Doppler Tissue Imaging Measurement and Calculation Labels

All DTI (Doppler Tissue Imaging) values in the patient report include prefix labels, **m** for medial or **l** for lateral, to indicate the specific mitral valve.

Measurement Label	Description
Med E' Vel	Early diastolic motion
AR/DR(m)	Begins the workflow for completing the following measurements: ▪ Med Accel Time - Acceleration Time of Med E' Vel ▪ Med Accel Rate - Acceleration Rate of Med E' Vel ▪ Med Decel Time - Deceleration Time of Med E' Vel ▪ Med Decel Rate - Deceleration Rate of Med E' Vel
Med A' Vel	Late diastolic motion
Med S Vel	Systolic motion
Med Accel Time	Med Acceleration Time
Med Decel Time	Med Deceleration Time
Lat E' Vel	Early diastolic motion
AR/DR(l)	Begins the workflow for completing the following measurements: ▪ Lat Accel Time - Acceleration Time of Lat E' Vel ▪ Lat Accel Rate - Acceleration Rate of Lat E' Vel ▪ Lat Decel Time - Deceleration Time of Lat E' Vel ▪ Lat Decel Rate - Deceleration Rate of Lat E' Vel
Lat A' Vel	Late diastolic motion
Lat S Vel	Systolic motion

Note: The system generates the calculation for E'/A'(Med) or E'/A'(Lat) when both Med E' Vel and Med A' Vel measurements are completed.

Sections in the Patient Report

Section Name	Description
Exam Specific Information	Indications and additional information for specific exam types
Assessment	Descriptions for exam-related structures. You can select system defined or user-defined descriptions
Stress Echo LV Results	(Available only for the cardiac exam type) Descriptions for stress echo measurement results for the left ventricle
WMS Results	(Available only for the cardiac exam type) Descriptions for wall motion scoring (WMS)
WMS Summary	(Available only for the cardiac exam type) Summary for wall motion scoring (WMS)
Images	Images included from the thumbnail panel
Comments	Any combination of predefined comments and user-entered comments
Measurement result by package	▪ Teichholz ▪ Cubed ▪ Gibson ▪ Simpson BP ▪ Simpson SP ▪ Mod. Simpson ▪ Single Plane ▪ Bi-Plane ▪ Bullet ▪ Mitral Valve ▪ DTI (MV) ▪ DTI (RV) ▪ MR ▪ PISA (AR) ▪ PISA (TR) ▪ PISA (MR) ▪ PISA (MS) ▪ AO/LA ▪ Aortic Valve ▪ AR ▪ Pulmonic Valve ▪ Pulmonary Vein ▪ LV Mass A-L ▪ LV Mass T-E ▪ LA Vol (A-L) ▪ LA Vol (Simp) ▪ RA Vol (A-L) ▪ RA Vol (Simp) ▪ RV ▪ Teichholz (M) ▪ Cubed (M) ▪ Gibson (M) ▪ AO/LA (M) ▪ Mitral Valve (M) ▪ RV (M) ▪ PR ▪ Tricuspid Valve ▪ TR ▪ Tei-index ▪ Qp/Qs

Appendix E Exam-specific Reference

Abdominal Measurements and Calculations.....	3
2D-Mode Measurement Labels	3
2D-Mode Calculation Labels.....	7
Spleen Volume Reference.....	8
Residual Urine Reference	8
M-Mode Measurement Labels	8
Doppler Measurement Labels.....	8
Doppler Calculation Labels.....	10
Sections in the Patient Report	11
Gynecology Measurements and Calculations.....	11
Gynecology 2D-Mode Measurement Labels.....	11
2D-Mode Calculation Labels.....	14
M-Mode Measurement Labels	15
Doppler Measurement Labels.....	15
Sections in the Patient Report	15
Pediatric Measurements and Calculations	15
2D-Mode Measurement Labels	15
2D-Mode Calculation Labels.....	19
Doppler Measurement Labels.....	19
Hip Angle Reference	21
Sections in the Patient Report	22
Neonatal Measurements and Calculations	22
Sections in the Patient Report	22
Vascular Measurements and Calculations	23
2D-Mode Measurement Labels	23
M-Mode Measurement Labels	26
Doppler Measurement Labels.....	26
Doppler Calculation Labels.....	29
Sections in the Patient Report	29

Urology Measurements and Calculations	30
2D-Mode Measurement Labels	30
2D-Mode Calculation Labels.....	31
Doppler Measurement Labels.....	31
Prostate Volume Reference	31
Prostate Specific Antigen Density Reference	31
Sections in the Patient Report	32
Musculoskeletal Measurements and Calculations	32
Sections in the Patient Report	32
Intraoperative Measurements and Calculations	33
Sections in the Patient Report	33
Small Parts Measurements and Calculations.....	33
2D-Mode Measurement Labels	33
2D-Mode Calculation Labels.....	34
Thyroid Lobe Volume Reference.....	34
Doppler Measurement Labels.....	35
M-Mode Calculation Labels	36
Sections in the Patient Report	36

Abdominal Measurements and Calculations

2D-Mode Measurement Labels

Measurement Label	Description	Measurement Method
Kidney Length	Kidney Length	Distance
Kidney Width	Kidney Width	Distance
Kidney AP	Kidney Anteroposterior	Distance
Cortical Thickness	Kidney Cortical Thickness	Distance
Liver Length	Liver Length	Distance
Liver AP	Liver Anteroposterior	Distance
Liver Width	Liver Width	Distance
Liver Mass 1	Liver Mass, numbered 1 through 10	3 Dist
CHD	Common Hepatic Duct	Distance
CBD	Common Bile Duct	Distance
Stent Portal End	Stent-Portal End (Proximal)	Distance
Stent Mid	Stent-Mid (Mid)	Distance
Stent Hepatic End	Stent-Hepatic End (Distal)	Distance
Donor A	Donor Artery	Distance
Portal V Pre Anast	Portal Vein Pre-anastomosis	Distance
Portal V Anast	Portal Vein Anastomosis	Distance
Transplant Kidney AP	Transplant Kidney Anteroposterior	Distance
Transplant Kidney Length	Transplant Kidney Length	Distance
Transplant Kidney Width	Transplant Kidney Width	Distance
GB AP	Gallbladder Anteroposterior	Distance
GB Length	Gallbladder Length	Distance
GB Width	Gallbladder Width	Distance
GB Wall Thickness	Gallbladder Wall Thickness	Distance
Pancreas Duct	Pancreas Duct	Distance
Pancreas Body	Pancreas Body	Distance
Pancreas Head	Pancreas Head	Distance
Pancreas Tail	Pancreas Tail	Distance
Spleen Length	Spleen Length	Distance
Spleen AP	Spleen Anteroposterior	Distance
Spleen Width	Spleen Width	Distance
Splenic A	Splenic Artery	Distance

Measurement Label	Description	Measurement Method
CIA	Common Iliac Artery	Distance
CIA AP	Common Iliac Artery Anteroposterior Diameter	Distance
CIA Trans	Common Iliac Artery Transverse Diameter	Distance
EIA	External Iliac Artery	Distance
Inflow Vessel 1	Inflow Vessel, numbered 1 through 2	Distance
Anast 1	Anastomosis, numbered 1 through 2	Distance
Graft 1	Graft, numbered 1 through 2	Distance
Outflow Vessel 1	Outflow Vessel, numbered 1 through 2	Distance
Art Anast	Art Anastomosis	Distance
EIV	External Iliac Vein	Distance
IIA	Internal Iliac Artery	Distance
IIA Trans	Internal Iliac Artery Transverse	Distance
IIV	Internal Iliac Vein	Distance
Aorta	Aorta	Distance
Aorta AP	Aorta Anteroposterior	Distance
Aorta Trans	Aorta Transverse	Distance
Suprarenal Ao	Suprarenal Aorta	Distance
Suprarenal Ao AP	Suprarenal Aorta Anteroposterior	Distance
Suprarenal Ao Diam	Suprarenal Aorta Diam	Distance
AAA	Abdominal Aortic Aneurysm	Distance
AAA AP	Abdominal Aortic Aneurysm Anteroposterior	Distance
AAA Trans	Abdominal Aortic Aneurysm Transverse	Distance
AAA dist from SMA	Abdominal Aortic Aneurysm distance from superior mesenteric artery	Distance
Neck Length	Infra-renal Neck Length	Distance
Neck Diam	Infra-renal Neck Diam	Distance
AAA Sac Length	Abdominal Aortic Aneurysm Sac Length	Distance
Ao Bif	Aorta at Bifurcation	Distance
Ao Bif AP	Aorta at Bifurcation Anteroposterior	Distance
Ao Bif Trans	Aorta at Bifurcation Transverse	Distance
CIA Stent	Common Iliac Artery Stent	Distance
CIA Stent AP	Common Iliac Artery Stent Anteroposterior	Distance
CIA Stent Trans	Common Iliac Artery Stent Transverse	Distance
CIA Aneurysm	Common Iliac Artery Aneurysm	Distance
CIA Aneurysm AP	Common Iliac Artery Aneurysm Anteroposterior	Distance
CIA Aneurysm Trans	Common Iliac Artery Aneurysm Transverse	Distance

Measurement Label	Description	Measurement Method
EIA Stent	External Iliac Artery Stent	Distance
EIA Stent AP	External Iliac Artery Stent Anteroposterior	Distance
EIA Stent Trans	External Iliac Artery Stent Transverse	Distance
IIA Stent	Internal Iliac Artery Stent	Distance
IIA Stent AP	Internal Iliac Artery Stent Anteroposterior	Distance
IIA Stent Trans	Internal Iliac Artery Stent Transverse	Distance
Prox Landing Site	Proximal (Stent Graft) Landing Site	Distance
Prox Landing Site AP	Proximal (Stent Graft) Landing Site Anteroposterior Diameter	Distance
Prox Landing Site Trans	Proximal (Stent Graft) Landing Site Transverse Diameter	Distance
Ao Stent Graft	Aortic Stent Graft	Distance
Ao Stent Graft AP	Aortic Stent Graft Anteroposterior Diameter	Distance
Ao Stent Graft Trans	Aortic Stent Graft Trans Diameter	Distance
Dist Landing Site	Distal (Stent Graft) Landing Site	Distance
Dist Landing Site AP	Distal (Stent Graft) Landing Site Anteroposterior Diameter	Distance
Dist Landing Site Trans	Distal (Stent Graft) Landing Site Transverse Diameter	Distance
Main RA Origin	Main Renal Artery Origin	Distance
IVC	Inferior Vena Cava	Distance
CIV	Common Iliac Vein	Distance
Celiac A Origin	Celiac Artery	Distance
SMA Origin	Superior Mesenteric Artery	Distance
SMA	Superior Mesenteric Artery	Distance
IMA Origin	Inferior Mesenteric Artery	Distance
IMA	Inferior Mesenteric Artery	Distance
Celiac A Origin, Postprandial	Celiac Artery-Postprandial	Distance
SMA Origin, Postprandial	Superior Mesenteric Artery-Postprandial	Distance
SMA Postprandial	Superior Mesenteric Artery-Postprandial	Distance
MPV	Main Portal Vein	Distance
LPV	Left Main Branch of Portal Vein	Distance
RPV	Right Main Branch of Portal Vein	Distance
SMV	Superior Mesenteric Vein	Distance
Splenic V	Splenic Vein	Distance

Measurement Label	Description	Measurement Method
CHA	Common Hepatic Artery	Distance
LHA	Left Branch of Hepatic Artery	Distance
RHA	Right Branch of Hepatic Artery	Distance
LHV	Left Hepatic Vein	Distance
MHV	Middle Hepatic Vein	Distance
RHV	Right Hepatic Vein	Distance
SVC	Superior Vena Cava	Distance
RA Stent	Renal Artery Stent	Distance
Main RA	Main Renal Artery	Distance
ARA Origin	Accessory Renal Artery Origin	Distance
ARA	Accessory Renal Artery	Distance
Renal V	Renal Vein	Distance
Transplant RA	Transplant Renal Artery	Distance
Transplant Renal V	Transplant Renal Vein	Distance
Renal V Anast	Renal Vein Anastomosis	Distance
Ureter	Ureter Length	Distance
Urethra	Urethra	Distance
PreVoid Bladder Length	Pre-Void Bladder Length	Distance
PreVoid Bladder AP	Pre-Void Bladder Anteroposterior	Distance
PreVoid Bladder Width	Pre-Void Bladder Width	Distance
PostVoid Bladder Length	Post-Void Bladder Length	Distance
PostVoid Bladder AP	Post-Void Bladder Anteroposterior	Distance
PostVoid Bladder Width	Post-Void Bladder Width	Distance
Bladder Length	Bladder Length	Distance
Bladder Width	Bladder Width	Distance
Bladder AP	Bladder Anteroposterior	Distance
Liver Segment 1*	Liver segments, numbered 1 through 8	---
Mass 1*	Mass, numbered 1 through 10	---
Site 1*	Site, numbered 1 through 10	---

*Available only for Virtual Touch applications.

2D-Mode Calculation Labels

Calculation Label	Description	Required Measurements	Units
Kidney Volume	Kidney Volume $3.1415926535 / 6 * \text{Kidney Length} * \text{Kidney Width} * \text{Kidney AP}$	Kidney Length Kidney Width Kidney AP	mL
Liver	Liver $3.1415926535 / 6 * \text{Liver Length} * \text{Liver Width} * \text{Liver AP}$	Liver Length Liver AP Liver Width	mL
Transplant Kidney Vol	Transplant Kidney Volume $3.1415926535 / 6 * \text{Transplant Kidney Length} * \text{Transplant Kidney Width} * \text{Transplant Kidney AP}$	Transplant Kidney Length Transplant Kidney AP Transplant Kidney Width	mL
GB	Gallbladder $3.1415926535 / 6 * \text{GB Length} * \text{GB Width} * \text{GB AP}$	GB Length GB Width GB AP	mL
Spleen Volume	Spleen Volume $3.1415926535 / 6 * \text{Spleen Length} * \text{Spleen AP} * \text{Spleen Width}$	Spleen Length Spleen AP Spleen Width	mL
EF	Ejection Fraction $\text{Stroke Volume} / \text{EDV} * 100$	EDV(Simpson SP) ESV(Simpson SP)	%
PreVoid Bladder Vol	PreVoid Bladder Volume $3.1415926535 / 6 * \text{PreVoid Bladder Length} * \text{PreVoid Bladder AP} * \text{PreVoid Bladder Width}$	PreVoid Bladder Length PreVoid Bladder AP PreVoid Bladder Width	mL
PostVoid Bladder Vol	PostVoid Bladder Volume $3.1415926535 / 6 * \text{PostVoid Bladder Length} * \text{PostVoid Bladder AP} * \text{PostVoid Bladder Width}$	PostVoid Bladder Length PostVoid Bladder AP PostVoid Bladder Width	mL
Micturated Volume	Micturated Volume $\text{PreVoid Bladder} - \text{PostVoid Bladder}$	PreVoid Bladder PostVoid Bladder	mL
Residual Urine (Haylen)	Residual Urine(Haylen) $\text{Bladder Width} * \text{Bladder AP} * 5.9 - 14.9$	Bladder Width Bladder AP	mL
Residual Urine (Dietz)	Residual Urine(Dietz) $\text{Bladder Length} * \text{Bladder AP} * 5.6$	Bladder Length Bladder AP	mL
Residual Urine (Dicuio)	Residual Urine(Dicuio) $\text{Bladder Length} * \text{Bladder AP} * \text{Bladder Width} * 0.5$	Bladder Length Bladder Width Bladder AP	mL

Spleen Volume Reference

Yetter EM, Acosta KB, Olson MC, Blundell K. Estimating Splenic Volume: "Sonographic Measurements Correlated with Helical CT Determination." *American Journal of Roentgenology* 181(6):1615–1620, 2003.

Residual Urine Reference

Haylen BT. "Verification of the Accuracy and Range of Transvaginal Ultrasound in Measuring Bladder Volumes in Women." *British Journal of Urology* 64(4):350–352, 1989.

Dietz HP, Velez D, Shek KL, Martin A. "Determination of Postvoid Residual by Translabial Ultrasound." *International Urogynecology Journal* 23(12):1749–1752, 2012.

Dicuio M, Pomara G, Menchini Fabris F, Ales V, Dahlstrand C, Morelli G. "Measurements of Urinary Bladder Volume: Comparison of Five Ultrasound Calculation Methods in Volunteers." *Archivio Italiano di Urologia e Andrologia* 77(1):60–62, 2005.

M-Mode Measurement Labels

Measurement Label	Description	Measurement Method
Fetal HR	Fetal Heart Rate, in beats per minute	HR

Doppler Measurement Labels

Measurement Label	Description	Measurement Method
CIA	Common Iliac Artery	PI Auto
EIA	External Iliac Artery	PI Auto
Inflow Vessel 1	Inflow Vessel, numbered 1 through 2	Velocity
Anast 1	Anastomosis, numbered 1 through 2	Velocity
Graft 1	Graft, numbered 1 through 2	Velocity
Outflow Vessel 1	Outflow Vessel, numbered 1 through 2	Velocity
Art Anast	Arterial Anastomosis	Velocity
EIV	External Iliac Vein	Velocity
IIA	Internal Iliac Artery	PI Auto
IIV	Internal Iliac Vein	Velocity
Aorta	Aorta	PI Auto
Suprarenal Ao	Suprarenal Aorta	PI Auto
AAA	Abdominal Aortic Aneurysm	RI
Ao Bif	Aorta at Bifurcation	PI Auto
CIA Stent	Common Iliac Artery Stent	PI Auto
CIA Aneurysm	Common Iliac Artery Aneurysm	RI
EIA Stent	External Iliac Artery Stent	RI
IIA Stent	Internal Iliac Artery Stent	RI

Measurement Label	Description	Measurement Method
Prox Landing Site	Proximal (Stent Graft) Landing Site	RI
Ao Stent Graft	Aortic Stent Graft	RI
Dist Landing Site	Distal (Stent Graft) Landing Site	RI
IVC	Inferior Vena Cava	Velocity
CIV	Common Iliac Vein	Velocity
Celiac A Origin	Celiac Artery Origin	PI Auto
SMA Origin	Superior Mesenteric Artery Origin	RI
SMA	Superior Mesenteric Artery	PI Auto
IMA Origin	Inferior Mesenteric Artery Origin	RI
IMA	Inferior Mesenteric Artery	PI Auto
Celiac A Origin, Postprandial	Celiac Artery-Postprandial	RI
SMA Origin, Postprandial	Superior Mesenteric Artery-Postprandial	RI
SMA Postprandial	Superior Mesenteric Artery-Postprandial	RI
MPV	Main Portal Vein	Velocity
LPV	Left Main Branch of Portal Vein	Velocity
RPV	Right Main Branch of Portal Vein	Velocity
SMV	Superior Mesenteric Vein	Velocity
Splenic V	Splenic Vein	Velocity
CHA	Common Hepatic Artery	RI
LHA	Left Branch of Hepatic Artery	RI
RHA	Right Branch of Hepatic Artery	RI
LHV	Left Hepatic Vein	Velocity
MHV	Middle Hepatic Vein	Velocity
RHV	Right Hepatic Vein	Velocity
Stent Portal End	Stent-Portal End (Proximal)	Velocity
Stent Mid	Stent-Mid (Mid)	Velocity
Stent Hepatic End	Stent-Hepatic End (Distal)	Velocity
Donor A	Donor Artery	RI
Portal V Pre Anast	Portal Vein Pre-anastomosis	Velocity
Portal V Anast	Portal Vein Anastomosis	Velocity
SVC	Superior Vena Cava	Velocity
RA Stent	Renal Artery Stent	RI
Main RA	Main Renal Artery	PI Auto
Main RA Origin	Main Renal Artery Origin	PI Auto
ARA Origin	Accessory Renal Artery Origin	RI
ARA	Accessory Renal Artery	RI

Measurement Label	Description	Measurement Method
Segmental A Sup	Superior Segmental Artery	RI
Segmental A Mid	Medial Segmental Artery	RI
Segmental A Inf	Inferior Segmental Artery	RI
Arcuate A Sup	Superior Arcuate Artery	RI
Arcuate A Mid	Medial Arcuate Artery	RI
Arcuate A Inf	Inferior Arcuate Artery	RI
Interlobar A Sup	Superior Interlobar Artery	RI
Interlobar A Mid	Medial Interlobar Artery	RI
Interlobar A Inf	Inferior Interlobar Artery	RI
Renal V	Renal Vein	Velocity
Transplant RA	Transplant Renal Artery	RI
Transplant Renal V	Transplant Renal Vein	Velocity
Renal V Anast	Renal Vein Anastomosis	Velocity
Splenic A	Splenic Artery	PI Auto
Fetal HR	Fetal Heart Rate, in beats per minute	HR

Doppler Calculation Labels

Calculation Label	Description	Required Measurements
Portal V Anast/Pre Anast Vel	Portal Vein Anastomosis / Portal Vein Anastomosis Velocity Ratio	Portal V Anast Pre Anast Vel
Main RA / Suprarenal Ao	Main Renal Artery / Suprarenal Aorta Ratio max(Main RA PS,Main RA PS Prox,Main RA PS Mid,Main RA PS Dist) / Suprarenal Aorta PS	Main RA PS Suprarenal Ao PS
Main RA Origin / Suprarenal Ao	Main Renal Artery Origin / Suprarenal Aorta Ratio max(Main RA PS,Main RA PS Prox,Main RA PS Mid,Main RA PS Dist) / Suprarenal Aorta PS	Main RA Origin Suprarenal Ao PS
ARA / Suprarenal Ao	Accessory Renal Artery / Suprarenal Aorta Ratio max(ARA PS,ARA PS Prox,ARA PS Mid,ARA PS Dist) / Suprarenal Aorta PS	ARA PS Suprarenal Ao PS
ARA Origin / Suprarenal Ao	Accessory Renal Artery Origin / Suprarenal Aorta Ratio max(ARA PS,ARA PS Prox,ARA PS Mid,ARA PS Dist) / Suprarenal Aorta PS	Suprarenal Ao PS ARA Origin

Sections in the Patient Report

Section Name	Description
Exam Specific Information	Indications and additional information for specific exam types
2D Measurements	Measurement labels and results for 2D-mode measurements
2D Mode Side by Side	Measurement labels and results for 2D-mode measurements for both right and left sides
M Measurements	Measurement labels and results for M-mode measurements
D Measurements	Measurement labels and results for Doppler measurements
D Mode Side by Side	Measurement labels and results for Doppler measurements for both right and left sides
Assessment	Descriptions for exam-related structures. You can select system defined or user-defined descriptions.
Images	Images included from the thumbnail panel
Comments	Any combination of predefined comments and user-entered comments

Gynecology Measurements and Calculations

Gynecology 2D-Mode Measurement Labels

Measurement Label	Description	Measurement Method
Cervix (TA)	Cervix Length, transabdominal	Distance
Cervix (EV)	Cervix Length, endovaginal	Distance
Kidney Length	Kidney Length	Distance
Kidney Width	Kidney Width	Distance
Kidney AP	Kidney AP	Distance
Ureter	Ureter Length	Distance
PreVoid Bladder Length	Pre-Void Bladder Length	Distance
PreVoid Bladder AP	Pre-Void Bladder Anteroposterior	Distance
PreVoid Bladder Width	Pre-Void Bladder Width	Distance
PostVoid Bladder Length	Post-Void Bladder Length	Distance
PostVoid Bladder AP	Post-Void Bladder Anteroposterior	Distance
PostVoid Bladder Width	Post-Void Bladder Width	Distance
Fibroid 1	Fibroid, numbered 1 through 5	Distance
Mass 1	Mass, numbered 1 through 5	Distance
Ovary Width (TA)	Ovary Width, transabdominal	Distance
Ovary Length (TA)	Ovary Length, transabdominal	Distance
Ovary AP (TA)	Ovary AP, transabdominal	Distance
Ovary Width (EV)	Ovarian Transverse Diameter (Width), endovaginal	Distance
Ovary Length (EV)	Ovarian Length, endovaginal	Distance
Ovary AP (EV)	Ovary Anteroposterior Diameter, endovaginal	Distance

Measurement Label	Description	Measurement Method
Endometrium	Endometrium	Distance
Endo Thick Width (TA)	Endometrium Thickness Width, transabdominal	Distance
Endo Thick Length (TA)	Endometrium Thickness Length, transabdominal	Distance
Endo Thick AP (TA)	Endometrium Thickness Anteroposterior Diameter, transabdominal	Distance
Endo Thick Width (EV)	Endometrium Thickness Width, endovaginal	Distance
Endo Thick Length (EV)	Endometrium Thickness Length, endovaginal	Distance
Endo Thick AP (EV)	Endometrium Thickness Anteroposterior Diameter, endovaginal	Distance
Uterus Width (TA)	Uterus Width, transabdominal	Distance
Uterus Length (TA)	Uterus Length, transabdominal	Distance
Uterus AP (TA)	Uterus Anteroposterior Diameter, transabdominal	Distance
Uterus Width (EV)	Uterus Width, endovaginal	Distance
Uterus Length (EV)	Uterus Length, endovaginal	Distance
Uterus AP (EV)	Uterus Anteroposterior Diameter, endovaginal	Distance
Uterine Width	Uterine Width	Distance
Uterine Length	Uterine Length	Distance
Uterine AP	Uterine Anteroposterior Diameter	Distance
Follicle 1	Follicle, numbered 1 through 15	2 Distance Avg
Cyst 1	Cyst, numbered 1 through 5	Distance
Bladder neck-Symphysis	Bladder neck-Symphysis Distance	Distance
Cervix descent	Cervix descent	Distance
Rectal Ampulla descent	Rectal Ampulla descent	Distance
Levator Hiatus AP	Levator Hiatus Anteroposterior	Distance
Levator Hiatus Lateral	Levator Hiatus Lateral	Distance
Levator Hiatus Area	Levator Hiatus Area	Ellipse
Rectovesical Angle (Rest)	Rectovesical Angle at Rest	Cross Angle
Rectovesical Angle (Valsalva)	Rectovesical Angle at Valsalva	Cross Angle
Detrusor wall thick.	Detrusor wall thickness	Distance
Urethral Rotation	Urethral Rotation	Cross Angle
RVA	Retrovesical Angle	Cross Angle
PUV	Posterior Urethrovesical Angle	Cross Angle

Measurement Label	Description	Measurement Method
Uterine Desc max	Maximum distance from the inferior margin of pubis symphysis to the cervix	Distance
Re.amp.desc.max	Maximum distance from the inferior margin of pubis symphysis to the rectal ampulla	Distance
Rectocele Depth	Rectocele is measured as the maximal Depth of the protrusion beyond the expected margin of the normal anterior rectal wall.	Distance
BND Rest	The bladder neck descent is the difference between the horizontal measurement of the bladder neck position at rest and at maximum valsalva	Distance
BND Stress	The bladder neck descent is the difference between the horizontal measurement of the bladder neck position at rest and at maximum valsalva	Distance
Bladder desc.max	Maximum distance from the inferior margin of pubis symphysis to the bladder	Distance
Levator Hiatus Stress	Quantification of pelvic floor muscle activity. Reduction of the mid-sagittal diameter of the levator hiatus is measured at the level of minimal hiatal dimensions	Distance
Bladder Length	Bladder Length	Distance
Bladder Width	Bladder Width	Distance
Bladder AP	Bladder Anterior Posterior	Distance
Hiatus AP	Hiatus AP length	Distance
Hiatus Lateral	Hiatus Lateral Diameter	Distance
PRM inner perimeter	Puborectalis muscle inner perimeter	Spline
Hiatus inner area	Pelvic floor hiatus inner area	Spline
Hiatus outer area	Pelvic floor hiatus outer area	Spline
Levator ani Angle	Levator ani angle	Cross Angle
IAS diameter	Internal Anal Sphincter Diameter	Distance
EAS diameter	External Anal Sphincter Diameter	Distance
Cervix Len	Cervical Length	Distance
CRL	Crown Rump Length	Distance
Yolk Sac	Yolk Sac Diameter	Distance
Mean Sac Diam	Mean Gestational Sac Diameter	3 Distance Avg
GSD	Gestational Sac Maximum Length	Distance

2D-Mode Calculation Labels

Calculation Label	Description	Required Measurements	Units
Kidney Volume	Kidney Volume $3.1415926535 / 6 * \text{Kidney Length} * \text{Kidney Width} * \text{Kidney AP}$	Kidney Length Kidney Width Kidney AP	mL
PreVoid Bladder Vol	PreVoid Bladder Volume $3.1415926535 / 6 * \text{PreVoid Bladder Length} * \text{PreVoid Bladder AP} * \text{PreVoid Bladder Width}$	PreVoid Bladder Length PreVoid Bladder AP PreVoid Bladder Width	mL
PostVoid Bladder Vol	PostVoid Bladder Volume $3.1415926535 / 6 * \text{PostVoid Bladder Length} * \text{PostVoid Bladder AP} * \text{PostVoid Bladder Width}$	PostVoid Bladder Length PostVoid Bladder AP PostVoid Bladder Width	mL
Micturated Volume	Micturated Volume PreVoid Bladder - PostVoid Bladder	PreVoid Bladder PostVoid Bladder	mL
Residual Urine (Haylen)	Residual Urine (Haylen) $\text{Bladder Width} * \text{Bladder AP} * 5.9 - 14.9$	Bladder Width Bladder AP	mL
Residual Urine (Dietz)	Residual Urine (Dietz) $\text{Bladder Length} * \text{Bladder AP} * 5.6$	Bladder Length Bladder AP	mL
Residual Urine (Dicuio)	Residual Urine (Dicuio) $\text{Bladder Length} * \text{Bladder AP} * \text{Bladder Width} * 0.5$	Bladder Length Bladder Width Bladder AP	mL
Ovary Volume (TA)	Ovary Volume (Transabdominal) $3.1415926535 / 6 * \text{Ovary Length (TA)} * \text{Ovary Width (TA)} * \text{Ovary AP (TA)}$	Ovary Length (TA) Ovary Width(TA) Ovary AP (TA)	mL
Ovary Volume (EV)	Ovary Volume (Endovaginal) $3.1415926535 / 6 * \text{Ovary Length (EV)} * \text{Ovary Width (EV)} * \text{Ovary AP (EV)}$	Ovary Length (EV) Ovary Width(EV) Ovary AP (EV)	mL
Uterine Volume	Uterine Volume $3.1415926535 / 6 * \text{Uterus Length} * \text{Uterus Width} * \text{Uterus AP}$	Uterine Length Uterine Width Uterine AP	mL
Uterus Volume (TA)	Uterus Volume (Transabdominal) $3.1415926535 / 6 * \text{Uterus Length(TA)} * \text{Uterus Width(TA)} * \text{Uterus AP(TA)}$	Uterus (TA) Length Uterus (TA) Width Uterus (TA) AP	mL
Uterus Volume (EV)	Uterus Volume (Endovaginal) $3.1415926535 / 6 * \text{Uterus Length(EV)} * \text{Uterus Width(EV)} * \text{Uterus AP(EV)}$	Uterus (EV) Length Uterus (EV) Width Uterus (EV) AP	mL
Endo Thick (TA)	Endometrium Thickness (Transabdominal) $3.1415926535 / 6 * \text{Endo Thick AP (TA)} * \text{Endo Thick Length (TA)} * \text{Endo Thick Width (TA)}$	Endo Thick AP (TA) Endo Thick Length (TA) Endo Thick Width (TA)	mL
Endo Thick (EV)	Endometrium Thickness (Endovaginal) $3.1415926535 / 6 * \text{Endo Thick AP (EV)} * \text{Endo Thick Length (EV)} * \text{Endo Thick Width (EV)}$	Endo Thick AP (EV) Endo Thick Length (EV) Endo Thick Width (EV)	mL

M-Mode Measurement Labels

Measurement Label	Description	Measurement Method
Fetal HR	Fetal Heart Rate, in beats per minute	HR

Doppler Measurement Labels

Measurement Label	Description	Measurement Method
Uterine A	Uterine Artery	PI Auto
Ovarian A	Ovarian Artery	PI Auto
Fetal HR	Fetal Heart Rate, in beats per minute	HR
Ovarian V	Ovarian Vein	Velocity

Sections in the Patient Report

Section Name	Description
Exam Specific Information	Indications and additional information for specific exam types
2D Measurements	Measurement labels and results for 2D-mode measurements
2D Mode Side by Side	Measurement labels and results for 2D-mode measurements for both right and left sides
D Measurements	Measurement labels and results for Doppler measurements
D Mode Side by Side	Measurement labels and results for Doppler measurements for both right and left sides
Assessment	Use the configuration settings to create user-defined anatomy structures and descriptions.
Images	Images included from the thumbnail panel
Comments	Any combination of predefined comments and user-entered comments

Pediatric Measurements and Calculations

2D-Mode Measurement Labels

Measurement Label	Description	Measurement Method
Kidney Length	Kidney Length	Distance
Kidney Width	Kidney Width	Distance
Kidney AP	Kidney AP	Distance
Cortical Thickness	Kidney Cortical Thickness	Distance
CIA	Common Iliac Artery	Distance
EIA	External Iliac Artery	Distance
Inflow Vessel 1	Inflow Vessel, numbered 1 through 2	Distance
Anast 1	Anastomosis, numbered 1 through 2	Distance
Graft 1	Graft, numbered 1 through 2	Distance
Outflow Vessel 1	Outflow Vessel, numbered 1 through 2	Distance

Measurement Label	Description	Measurement Method
Art Anast	Art Anastomosis	Distance
EIV	External Iliac Vein	Distance
IIA	Internal Iliac Artery	Distance
IIV	Internal Iliac Vein	Distance
Aorta	Aorta Diameter	Distance
Suprarenal Ao	Suprarenal Aorta	Distance
Suprarenal Ao AP	Suprarenal Aorta AP	Distance
Suprarenal Ao Diam	Suprarenal Aorta Diam	Distance
AAA	Abdominal Aortic Aneurysm	Distance
AAA AP	Abdominal Aortic Aneurysm AP	Distance
AAA Trans	Abdominal Aortic Aneurysm Trans	Distance
AAA dist from SMA	Abdominal Aortic Aneurysm distance from SMA	Distance
Neck Length	Infra-renal Neck Length	Distance
Neck Diam	Infra-renal Neck Diam	Distance
AAA Sac Length	Abdominal Aortic Aneurysm Sac Length	Distance
Ao Bif	Aorta at Bifurcation	Distance
Ao Bif AP	Aorta at Bifurcation	Distance
Ao Bif Trans	Aorta at Bifurcation	Distance
CIA Stent	Common Iliac Artery Stent	Distance
CIA Stent AP	Common Iliac Artery Stent AP	Distance
CIA Stent Trans	Common Iliac Artery Stent Trans	Distance
CIA Aneurysm	Common Iliac Artery Aneurysm	Distance
CIA Aneurysm AP	Common Iliac Artery Aneurysm AP	Distance
CIA Aneurysm Trans	Common Iliac Artery Aneurysm Trans	Distance
EIA Stent	External Iliac Artery Stent	Distance
EIA Stent AP	External Iliac Artery Stent	Distance
EIA Stent Trans	External Iliac Artery Stent	Distance
IIA Stent	Internal Iliac Artery Stent	Distance
IIA Stent AP	Internal Iliac Artery Stent AP	Distance
IIA Stent Trans	Internal Iliac Artery Stent Trans	Distance
Prox Landing Site	Proximal (Stent Graft) Landing Site	Distance
Prox Landing Site AP	Proximal (Stent Graft) Landing Site AP Diameter	Distance
Prox Landing Site Trans	Proximal (Stent Graft) Landing Site Trans Diameter	Distance
Ao Stent Graft	Aortic Stent Graft	Distance
Ao Stent Graft AP	Aortic Stent Graft AP Diameter	Distance
Ao Stent Graft Trans	Aortic Stent Graft Trans Diameter	Distance

Measurement Label	Description	Measurement Method
Dist Landing Site	Distal (Stent Graft) Landing Site	Distance
Dist Landing Site AP	Distal (Stent Graft) Landing Site AP Diameter	Distance
Dist Landing Site Trans	Distal (Stent Graft) Landing Site Trans Diameter	Distance
Main RA Origin	Main Renal Artery Origin	Distance
IVC	Inferior Vena Cava	Distance
CIV	Common Iliac Vein	Distance
Celiac A Origin	Celiac Artery	Distance
SMA Origin	Superior Mesenteric Artery	Distance
SMA	Superior Mesenteric Artery	Distance
IMA Origin	Inferior Mesenteric Artery	Distance
IMA	Inferior Mesenteric Artery	Distance
Celiac A Origin, Postprandial	Celiac Artery-Postprandial	Distance
SMA Origin, Postprandial	Superior Mesenteric Artery-Postprandial	Distance
SMA Postprandial	Superior Mesenteric Artery-Postprandial	Distance
MPV	Main Portal Vein	Distance
LPV	Left Main Branch of Portal Vein	Distance
RPV	Right Main Branch of Portal Vein	Distance
SMV	Superior Mesenteric Vein	Distance
Splenic V	Splenic Vein	Distance
CHA	Common Hepatic Artery	Distance
LHA	Left Branch of Hepatic Artery	Distance
RHA	Right Branch of Hepatic Artery	Distance
LHV	Left Hepatic Vein	Distance
MHV	Middle Hepatic Vein	Distance
RHV	Right Hepatic Vein	Distance
Liver Length	Liver Length	Distance
Liver AP	Liver AP	Distance
Liver Width	Liver Width	Distance
Liver Mass 1	Liver Mass, numbered 1 through 10	3 Dist
CHD	Common Hepatic Duct	Distance
CBD	Common Bile Duct	Distance
Stent Portal End	Stent-Portal End (Proximal)	Distance
Stent Mid	Stent-Mid (Mid)	Distance
Stent Hepatic End	Stent-Hepatic End (Distal)	Distance

Measurement Label	Description	Measurement Method
Donor A	Donor Artery	Distance
Portal V Pre Anast	Portal Vein Pre-anastomosis	Distance
Portal V Anast	Portal Vein Anastomosis	Distance
SVC	Superior Vena Cava	Distance
RA Stent	Renal Artery Stent	Distance
Main RA	Main Renal Artery	Distance
ARA Origin	Accessory Renal Artery Origin	Distance
ARA	Accessory Renal Artery	Distance
Renal V	Renal Vein	Distance
Art Anast	Arterial Anastomosis	Distance
Transplant RA	Transplant Renal Artery	Distance
Transplant Renal V	Transplant Renal Vein	Distance
Renal V Anast	Renal Vein Anastomosis	Distance
Transplant Kidney AP	Transplant Kidney AP	Distance
Transplant Kidney Length	Transplant Kidney Length	Distance
Transplant Kidney Width	Transplant Kidney Width	Distance
GB AP	Gallbladder AP	Distance
GB Length	Gallbladder Length	Distance
GB Width	Gallbladder Width	Distance
GB Wall Thickness	Gallbladder Wall Thickness	Distance
Pancreas Duct	Pancreas Duct	Distance
Pancreas Body	Pancreas Body	Distance
Pancreas Head	Pancreas Head	Distance
Pancreas Tail	Pancreas Tail	Distance
Spleen Length	Spleen Length	Distance
Spleen AP	Spleen AP	Distance
Spleen Width	Spleen Width	Distance
Splenic A	Splenic Artery	Distance
Ureter	Ureter Length	Distance
Hip Angle	Hip angle	Angle
Lt Hip Angle	Left hip angle	Angle
Rt Hip Angle	Right hip angle	Angle
Flexed FHC	Flexed percentage of femoral head coverage	Ratio
Lt Flexed FHC	Left flexed percentage of femoral head coverage	Ratio
Rt Flexed FHC	Right flexed percentage of femoral head coverage	Ratio
Stressed FHC	Stressed percentage of femoral head coverage	Ratio
Lt Stressed FHC	Left stressed percentage of femoral head coverage	Ratio
Rt Stressed FHC	Right stressed percentage of femoral head coverage	Ratio

2D-Mode Calculation Labels

Calculation Label	Description	Required Measurements	Units
Kidney Volume	Kidney Volume	Kidney Length Kidney Width Kidney AP	mL
Liver	Liver	Liver Length Liver AP Liver Width	mL
GB	Gallbladder	GB Length GB Width GB AP	mL

Doppler Measurement Labels

Measurement Label	Description	Measurement Method
CIA	Common Iliac Artery	PI Auto
EIA	External Iliac Artery	PI Auto
Inflow Vessel 1	Inflow Vessel, numbered 1 through 2	Velocity
Anast 1	Anastomosis, numbered 1 through 2	Velocity
Graft 1	Graft, numbered 1 through 2	Velocity
Outflow Vessel 1	Outflow Vessel, numbered 1 through 2	Velocity
Art Anast	Art Anastomosis	Velocity
EIV	External Iliac Vein	Velocity
IIA	Internal Iliac Artery	PI Auto
IIV	Internal Iliac Vein	Velocity
Aorta	Intimal Diameter	PI Auto
Suprarenal Ao	Suprarenal Aorta	PI Auto
AAA	Abdominal Aortic Aneurysm	RI
Ao Bif	Aorta at Bifurcation	PI Auto
CIA Stent	Common Iliac Artery Stent	PI Auto
CIA Aneurysm	Common Iliac Artery Aneurysm	RI
EIA Stent	External Iliac Artery Stent	RI
IIA Stent	Internal Iliac Artery Stent	RI
Prox Landing Site	Proximal (Stent Graft) Landing Site	RI
Ao Stent Graft	Aortic Stent Graft	RI
Dist Landing Site	Distal (Stent Graft) Landing Site	RI

Measurement Label	Description	Measurement Method
IVC	Inferior Vena Cava	Velocity
CIV	Common Iliac Vein	Velocity
Celiac A Origin	Celiac Artery	PI Auto
SMA Origin	Superior Mesenteric Artery	RI
SMA	Superior Mesenteric Artery	PI Auto
IMA Origin	Inferior Mesenteric Artery	RI
IMA	Inferior Mesenteric Artery	PI Auto
Celiac A Origin, Postprandial	Celiac Artery-Postprandial	RI
SMA Origin, Postprandial	Superior Mesenteric Artery-Postprandial	RI
SMA Postprandial	Superior Mesenteric Artery-Postprandial	RI
MPV	Main Portal Vein	Velocity
LPV	Left Main Branch of Portal Vein	Velocity
RPV	Right Main Branch of Portal Vein	Velocity
SMV	Superior Mesenteric Vein	Velocity
Splenic V	Splenic Vein	Velocity
CHA	Common Hepatic Artery	RI
LHA	Left Branch of Hepatic Artery	RI
RHA	Right Branch of Hepatic Artery	RI
LHV	Left Hepatic Vein	Velocity
MHV	Middle Hepatic Vein	Velocity
RHV	Right Hepatic Vein	Velocity
Stent Portal End	Stent-Portal End (Proximal)	Velocity
Stent Mid	Stent-Mid (Mid)	Velocity
Stent Hepatic End	Stent-Hepatic End (Distal)	Velocity
Donor A	Donor Artery	RI
Portal V Pre Anast	Portal Vein Pre-anastomosis	Velocity
Portal V Anast	Portal Vein Anastomosis	Velocity
SVC	Superior Vena Cava	Velocity
RA Stent	Renal Artery Stent	RI
Main RA	Main Renal Artery	PI Auto
ARA Origin	Accessory Renal Artery Origin	RI
ARA	Accessory Renal Artery	RI

Measurement Label	Description	Measurement Method
Segmental A Sup	Superior Segmental Artery	RI
Segmental A Mid	Medial Segmental Artery	RI
Segmental A Inf	Inferior Segmental Artery	RI
Arcuate A Sup	Superior Arcuate Artery	RI
Arcuate A Mid	Medial Arcuate Artery	RI
Arcuate A Inf	Inferior Arcuate Artery	RI
Renal V	Renal Vein	Velocity
Art Anast	Arterial Anastomosis	Velocity
Transplant RA	Transplant Renal Artery	RI
Transplant Renal V	Transplant Renal Vein	Velocity
Renal V Anast	Renal Vein Anastomosis	Velocity
Splenic A	Splenic Artery	PI Auto

Hip Angle Reference

- Castro RH, Castro M, Fernandes TC, Braga A, Madureira AJBS. "Infant hip ultrasound – all you wanted to know but were afraid to ask." *European Society of Radiology*, 1-20, 2011.
- Clarke NMP, Harcke HT, McHugh P, Lee MS, Borns PF, MacEwen GD. "Real-time Ultrasound in the Diagnosis of Congenital Dislocation and Dysplasia of the Hip." *The Journal of Bone and Joint Surgery* 67(3):406–412, 1985.
- Graf R. "Sonographic Diagnosis of Hip Dysplasia and Hip Dislocation." Revised from reprint: Schuler P, Graf R. "Sonographie in der Orthopädie" in Braun - Günther - Schwerek. *Ultraschalldiagnostik*, 4. Erg. Lfg 7, 1986.
- Grissom LE, Harcke HT, Kumar SJ, Bassett GS, MacEwen GD. "Ultrasound Evaluation of Hip Position in the Pavlik Harness." *The American Institute of Ultrasound in Medicine* 7(1):1–6, 1988.
- Harcke HT, Grissom LE. "Performing Dynamic Sonography of the Infant Hip." *American Journal of Roentgenology* 155(4):837–844, 1990.
- Morin C, Harcke HT, MacEwen GD. "The Infant Hip: Real-Time US Assessment of Acetabular Development." *Radiology* 157(3):673-677, 1985.
- Peterlein CD, Schüttler KF, Lakemeier S, Timmesfeld N, Görg C, Fuchs-Winkelmann S, Schofer MD. "Reproducibility of Different Screening Classifications in Ultrasonography of the Newborn Hip." *BMC Pediatrics* 10:98, 2010.
- Smergel E, Losik SB, Rosenberg HK. "Sonography of Hip Dysplasia." *Ultrasound Quarterly* 20(4):201–216, 2004.

Sections in the Patient Report

Section Name	Description
Exam Specific Information	Indications and additional information for specific exam types
2D Measurements	Measurement labels and results for 2D-mode measurements
2D Mode Side by Side	Measurement labels and results for 2D-mode measurements for both right and left sides
M Measurements	Measurement labels and results for M-mode measurements
D Measurements	Measurement labels and results for Doppler measurements
D Mode Side by Side	Measurement labels and results for Doppler measurements for both right and left sides
Assessment	Descriptions for exam-related structures. You can select system defined or user-defined descriptions.
Images	Images included from the thumbnail panel
Comments	Any combination of predefined comments and user-entered comments

Neonatal Measurements and Calculations

The neonatal exam type has only general measurements and no measurement labels.

Sections in the Patient Report

Section Name	Description
Exam Specific Information	Includes systole and diastole for both left and right sides and the blood pressure for each side
2D Measurements	Measurement labels and results for 2D-mode measurements
2D Mode Side by Side	Measurement labels and results for 2D-mode measurements for both right and left sides
M Measurements	Measurement labels and results for M-mode measurements
D Measurements	Measurement labels and results for Doppler measurements
D Mode Side by Side	Measurement labels and results for Doppler measurements for both right and left sides
Assessment	Use the system configuration settings to create user-defined anatomy structures and descriptions.
Images	Images included from the thumbnail panel
Comments	Any combination of predefined comments and user-entered comments

Vascular Measurements and Calculations

2D-Mode Measurement Labels

Measurement Label	Description	Measurement Method
EIV	External Iliac Vein	Distance
IIV	Internal Iliac Vein	Distance
CIV	Common Iliac Vein	Distance
Axillary Vein	Axillary Vein	Distance
Brachial Vein	Brachial Vein	Distance
CCA	Common Carotid Artery	ET-%Stenosis
ICA	Internal Carotid Artery	ET-%Stenosis
ECA	External Carotid Artery	ET-%Stenosis
Vert A	Vertebral Artery	ET-%Stenosis
Bulb	Carotid Bulb	Distance
Subclav A	Subclavian Artery	Distance
Innom A	Innominate Artery	Distance
Bulb IMT	Carotid Bulb Intima-Media Thickness	---
CCA IMT	Common Carotid Artery Intima-Media Thickness	---
ICA IMT	Internal Carotid Artery Intima-Media Thickness	---
HR	Heart Rate	(user-entered value)
CFV	Common Femoral Vein	Distance
Deep Fem V	Deep Femoral Vein	Distance
Fem V	Femoral Vein	Distance
Pop V	Popliteal Vein	Distance
Med Gastroc V 1	Medial Gastrocnemius Vein, numbered 1 through 2	Distance
Lat Gastroc V 1	Lateral Gastrocnemius Vein, numbered 1 through 2	Distance
Soleus V	Soleus Vein	Distance
PTV 1	Posterior Tibial Vein, numbered 1 through 2	Distance
Peroneal V 1	Peroneal Vein, numbered 1 through 2	Distance
ATV 1	Anterior Tibial Vein, numbered 1 through 2	Distance
SFJ	Saphenofemoral Junction	Distance
GSV Upper Thigh	Great Saphenous Vein: Upper Thigh	Distance
GSV Mid Thigh	Great Saphenous Vein: Mid Thigh	Distance
GSV Lower Thigh	Great Saphenous Vein: Lower Thigh	Distance
GSV Knee	Great Saphenous Vein: Knee	Distance
GSV Upper Calf	Great Saphenous Vein: Upper Calf	Distance
GSV Mid Calf	Great Saphenous Vein: Mid Calf	Distance
GSV Lower Calf	Great Saphenous Vein: Lower Calf	Distance
GSV Ankle	Great Saphenous Vein: Ankle	Distance

Measurement Label	Description	Measurement Method
SPJ	Saphenopopliteal Junction	Distance
SSV Upper Calf	Small Saphenous Vein: Upper Calf	Distance
SSV Mid Calf	Small Saphenous Vein: Mid Calf	Distance
SSV Lower Calf	Small Saphenous Vein: Lower Calf	Distance
SSV Ankle	Small Saphenous Vein: Ankle	Distance
Cockett's Perf 1	Cockett's Perforator, numbered 1 through 3	Distance
Boyd's Perf	Boyd's Perforator	Distance
Dodd Perf	Dodd Perforator	Distance
Hunterian Perf	Hunterian Perforator	Distance
Innom V	Innominate Vein	Distance
IJV	Internal Jugular Vein	Distance
SCL V	Subclavian Vein	Distance
Brachial V1	Brachial Vein, numbered 1 through 2	Distance
Cephalic V	Cephalic Vein	Distance
Basilic V	Basilic Vein	Distance
Mcub V	Median Cubital Vein	Distance
Radial V	Radial Vein	Distance
Radial V1	Radial Vein, numbered 1 through 2	Distance
Ulnar V	Ulnar Vein	Distance
Ulnar V1	Ulnar Vein, numbered 1 through 2	Distance
MCA	Middle Cerebral Artery	Distance
ICA-Siphon	Internal Carotid Artery-Siphon	Distance
ACA A1	The A-1 section of the Anterior Cerebral Artery	Distance
ACA A2	The A-2 section of the Anterior Cerebral Artery	Distance
ACoA	Anterior Communicating Artery	Distance
PCA-P1	The P1 section of the Posterior Cerebral Artery	Distance
PCA-P2	The P2 section of the Posterior Cerebral Artery	Distance
PCoA	Posterior Communicating Artery	Distance
PCA	Posterior Cerebral Artery	Distance
Basilar A	Basilar Artery	Distance
Ophthalmic A	Ophthalmic Artery	Distance
CIA	Common Iliac Artery	Distance
EIA	External Iliac Artery	Distance

Measurement Label	Description	Measurement Method
Inflow Vessel 1	Inflow Vessel, numbered 1 through 2	Distance
Anast 1	Anastomosis, numbered 1 through 2	Distance
Graft 1	Graft, numbered 1 through 2	Distance
Outflow Vessel 1	Outflow Vessel, numbered 1 through 2	Distance
Art Anast	Arterial Anastomosis	Distance
Axillary A	Axillary Artery	Distance
Brachial A	Brachial Artery	Distance
CFA	Common Femoral Artery	Distance
CFA Pseudo	Common Femoral Artery Pseudoaneurysm	Distance
CFA Stent	Common Femoral Artery Stent	Distance
DFA	Deep Femoral Artery	Distance
DFA Stent	Deep Femoral Artery Stent	Distance
SFA	Superficial Femoral Artery	Distance
SFA Stent	Superficial Femoral Artery Stent	Distance
Pop A	Popliteal Artery	Distance
Pop A Aneurysm	Popliteal Artery Aneurysm	Distance
Pop A Aneurysm AP	Popliteal Artery Aneurysm	Distance
Pop A Aneurysm Trans	Popliteal Artery Aneurysm	Distance
Pop A Stent	Popliteal Artery Stent	Distance
T-P Trunk	Tibioperoneal Trunk	Distance
ATA	Anterior Tibial Artery	Distance
PTA	Posterior Tibial Artery	Distance
Peroneal A	Peroneal Artery	Distance
DPA	Dorsalis Pedis Artery	Distance
UE stent	Upper Extremity Stent	Distance
Deep Brachial A	Deep Brachial Artery	Distance
Radial A	Radial Artery	Distance
Ulnar A	Ulnar Artery	Distance
Inflo A	Inflow Artery	Distance
AVF	Arteriovenous Fistula	Distance
Venous Anast	Venous Anastomosis	Distance
Outflow V	Outflow Vein	Distance
Inflow Vessel	Inflow Vessel	Distance
Anast	Anastomosis	Distance
Graft	Graft	Distance
Outflow Vessel	Outflow Vessel	Distance

M-Mode Measurement Labels

Measurement Label	Description	Measurement Method
HR	Heart Rate	HR

Doppler Measurement Labels

Measurement Label	Description	Measurement Method
IIV	Internal Iliac Vein	Velocity
CIV	Common Iliac Vein	Velocity
Axillary Vein	Axillary Vein	Velocity
Brachial Vein	Brachial Vein	Velocity
CCA	Common Carotid Artery	RI
ICA	Internal Carotid Artery	RI
ECA	External Carotid Artery	RI
Vert A	Vertebral Artery	RI
Bulb	Carotid Bulb	RI
Subclav A	Subclavian Artery	RI
Innom A	Innominate Artery	RI
HR	Heart Rate	HR
CFV	Common Femoral Vein	Velocity
Deep Fem V	Deep Femoral Vein	Velocity
Fem V	Femoral Vein	Velocity
Pop V	Popliteal Vein	Velocity
Med Gastroc V 1	Medial Gastrocnemius Vein, numbered 1 through 2	Velocity
Lat Gastroc V 1	Lateral Gastrocnemius Vein, numbered 1 through 2	Velocity
Soleus V	Soleus Vein	Velocity
PTV 1	Posterior Tibial Vein, numbered 1 through 2	Velocity
Peroneal V 1	Peroneal Vein, numbered 1 through 2	Velocity
ATV 1	Anterior Tibial Vein, numbered 1 through 2	Velocity
SFJ	Saphenofemoral Junction	Velocity
GSV Upper Thigh	Great Saphenous Vein: Upper Thigh	Velocity
GSV Mid Thigh	Great Saphenous Vein: Mid Thigh	Velocity
GSV Lower Thigh	Great Saphenous Vein: Lower Thigh	Velocity
GSV Knee	Great Saphenous Vein: Knee	Velocity
GSV Upper Calf	Great Saphenous Vein: Upper Calf	Velocity
GSV Mid Calf	Great Saphenous Vein: Mid Calf	Velocity
GSV Lower Calf	Great Saphenous Vein: Lower Calf	Velocity
GSV Ankle	Great Saphenous Vein: Ankle	Velocity

Measurement Label	Description	Measurement Method
SPJ	Saphenopopliteal Junction	Velocity
SSV Upper Calf	Small Saphenous Vein: Upper Calf	Velocity
SSV Mid Calf	Small Saphenous Vein: Mid Calf	Velocity
SSV Lower Calf	Small Saphenous Vein: Lower Calf	Velocity
SSV Ankle	Small Saphenous Vein: Ankle	Velocity
Cockett's Perf 1	Cockett's Perforator, numbered 1 through 3	Velocity
Boyd's Perf	Boyd's Perforator	Velocity
Dodd Perf	Dodd Perforator	Velocity
Hunterian Perf	Hunterian Perforator	Velocity
Innom V	Innominate Vein	Velocity
IJV	Internal Jugular Vein	Velocity
SCL V	Subclavian Vein	Velocity
Brachial V1	Brachial Vein, numbered 1 through 2	Velocity
Cephalic V	Cephalic Vein	Velocity
Basilic V	Basilic Vein	Velocity
Mcub. V	Median Cubital Vein	Velocity
Radial V	Radial Vein	Velocity
Radial V1	Radial Vein, numbered 1 through 2	Velocity
Ulnar V	Ulnar Vein	Velocity
Ulnar V1	Ulnar Vein, numbered 1 through 2	Velocity
Dist ICA skull entrance	Distal Internal Carotid Artery at entrance to skull	RI
Terminal ICA	Terminal Internal Carotid Artery	RI
Terminal Vert A	Terminal Vertebral Artery	RI
MCA	Middle Cerebral Artery	Velocity
ICA-Siphon	Internal Carotid Artery-Siphon	Velocity
ACA A1	The A-1 section of the Anterior Cerebral Artery	Velocity
ACA A2	The A-2 section of the Anterior Cerebral Artery	Velocity
ACoA	Anterior Communicating Artery	Velocity
PCA P1	The P1 section of the Posterior Cerebral Artery	Velocity
PCA P2	The P2 section of the Posterior Cerebral Artery	Velocity
PCoA	Posterior Communicating Artery	Velocity
PCA	Posterior Cerebral Artery	Velocity
Basilar A	Basilar Artery	Velocity
M1 MCA	Middle Cerebral Artery M1 Segment	Velocity
M2 MCA	Middle Cerebral Artery M2 Segment	Velocity
ACA	Anterior Cerebral Artery	Velocity

Measurement Label	Description	Measurement Method
EIV	External Iliac Vein	Velocity
CIA	Common Iliac Artery	PI Auto
EIA	External Iliac Artery	PI Auto
Inflow Vessel 1	Inflow Vessel, numbered 1 through 2	Velocity
Anast 1	Anastomosis, numbered 1 through 2	Velocity
Graft 1	Graft, numbered 1 through 2	Velocity
Outflow Vessel 1	Outflow Vessel, numbered 1 through 2	Velocity
Art Anast	Arterial Anastomosis	RI
Axillary A	Axillary Artery	PI Auto
Brachial A	Brachial Artery	PI Auto
Subclav A	Subclavian Artery	RI
Innom A	Innominate Artery	RI
CFA	Common Femoral Artery	Velocity
CFA Pseudo	Common Femoral Artery Pseudoaneurysm	Velocity
CFA Stent	Common Femoral Artery Stent	Velocity
DFA	Deep Femoral Artery	Velocity
DFA Stent	Deep Femoral Artery Stent	Velocity
SFA	Superficial Femoral Artery	Velocity
SFA Stent	Superficial Femoral Artery Stent	Velocity
Pop A	Popliteal Artery	Velocity
Pop A Aneurysm	Popliteal Artery Aneurysm	Velocity
Pop A Stent	Popliteal Artery Stent	Velocity
T-P Trunk	Tibioperoneal Trunk	Velocity
ATA	Anterior Tibial Artery	Velocity
PTA	Posterior Tibial Artery	Velocity
Peroneal A	Peroneal Artery	Velocity
DPA	Dorsalis Pedis Artery	Velocity
UE Stent	Upper Extremity Stent	Velocity

Measurement Label	Description	Measurement Method
Deep Brachial A	Deep Brachial Artery	Velocity
Radial A	Radial Artery	Velocity
Ulnar A	Ulnar Artery	Velocity
Inflow A	Inflow Artery	Velocity
AVF	Arteriovenous Fistula	Velocity
Venous Anast	Venous Anastomosis	Velocity
Outflow V	Outflow Vein	Velocity
Inflow Vessel	Inflow Vessel	Velocity
Anast	Anastomosis	Velocity
Outflow Vessel	Outflow Vessel	Velocity

Doppler Calculation Labels

Calculation Label	Description	Required Measurements
ICA / CCA (d)	Internal Carotid Artery / Common Carotid Artery (diastole)	ICA ED CCA ED
ICA / CCA (s)	Internal Carotid Artery / Common Carotid Artery (systole)	ICA PS CCA PS
MCA/ICA-Siphon	Middle Cerebral Artery / Internal Carotid Artery-Siphon	MCA ICA-Siphon

Sections in the Patient Report

Section Name	Description
Exam Specific Information	Indications and additional information for specific exam types
2D Measurements	Measurement labels and results for 2D-mode measurements
2D Mode Side by Side	Measurement labels and results for 2D-mode measurements for both right and left sides
IMT	(Available only for Intima Media Thickness measurements) Measurement labels and results for Intima Media Thickness (IMT)
IMT Side by Side	(Available only for Intima Media Thickness measurements) Measurement labels and results for Intima Media Thickness (IMT) for both right and left sides
M Measurements	Measurement labels and results for M-mode measurements
D Measurements	Measurement labels and results for Doppler measurements
D Mode Side by Side	Measurement labels and results for Doppler measurements for both right and left sides
Venous Assessment	Includes system-defined descriptions and observations for venous structures.
Carotid Assessment	Includes system-defined descriptions and observations for carotid structures.
Arterial Ratios	Includes arterial ratios
Images	Images included from the thumbnail panel
Comments	Any combination of predefined comments and user-entered comments

Urology Measurements and Calculations

2D-Mode Measurement Labels

Measurement Label	Description	Measurement Method
IIA	Internal Iliac Artery	Distance
IIV	Internal Iliac Vein	Distance
Ureter	Ureter Length	Distance
Urethra	Urethra	Distance
PreVoid Bladder Length	Pre-Void Bladder Length	Distance
PreVoid Bladder AP	Pre-Void Bladder Anteroposterior	Distance
PreVoid Bladder Width	Pre-Void Bladder Width	Distance
PostVoid Bladder Length	Post-Void Bladder Length	Distance
PostVoid Bladder AP	Post-Void Bladder Anteroposterior	Distance
PostVoid Bladder Width	Post-Void Bladder Width	Distance
Prostate V1	Prostate V1, numbered 1 through 2	3Dist Volume
Prostate V1 Trans	Prostate V1 Transverse, numbered 1 through 2	Distance
Prostate V1 Length	Prostate V1 Length, numbered 1 through 2	Distance
Prostate V1 AP	Prostate V1 Depth, numbered 1 through 2	Distance
Seminal Ves	Seminal Vesicle	Distance
Vas Deferens	Vas Deferens	Distance
Prostatic Urethra	Prostatic Urethra	Distance
Ejaculatory Duct	Ejaculatory Duct	Distance
Int Pudendal A	Internal Pudendal Artery	Distance
Pudendal V	Pudendal Vein	Distance
Mass 1	Urology Mass, numbered 1 through 3	Distance

2D-Mode Calculation Labels

Calculation Label	Description	Required Measurements	Units
PreVoid Bladder Vol	PreVoid Bladder Volume $3.1415926535 / 6 * \text{PreVoid Bladder Length} * \text{PreVoid Bladder AP} * \text{PreVoid Bladder Width}$	PreVoid Bladder Length PreVoid Bladder AP PreVoid Bladder Width	mL
PostVoid Bladder Vol	PostVoid Bladder Volume $3.1415926535 / 6 * \text{PostVoid Bladder Length} * \text{PostVoid Bladder AP} * \text{PostVoid Bladder Width}$	PostVoid Bladder Length PostVoid Bladder AP PostVoid Bladder Width	mL
Micturated Volume	Micturated Volume PreVoid Bladder - PostVoid Bladder	PreVoid Bladder PostVoid Bladder	mL
Prostate V1	Prostate Volume1	Prostate V1 Length Prostate V1 Width Prostate V1 AP	mL
Prostate V2	Prostate Volume2	Prostate V2 Length Prostate V2 Width Prostate V2 AP	mL
PSAD 1	Prostate Specific Antigen Density 1	Prostate V1 PSA	-
PSAD 2	Prostate Specific Antigen Density 2	Prostate V2 PSA	-

Doppler Measurement Labels

Measurement Label	Description	Measurement Method
IIA	Internal Iliac Artery	PI Auto
IIV	Internal Iliac Vein	Velocity
Int Pudendal A	Internal Pudendal Artery	Velocity
Pudendal V	Pudendal Vein	Velocity

Prostate Volume Reference

Rifkin MD. "Prostrate and Seminal Vesicle Measurements." Chapter 16 in Atlas of Ultrasound Measurements. B.B. Goldberg and A.B. Kurtz Chicago: Year Book Medical Publisher. 1990.

Prostate Specific Antigen Density Reference

Benson MC, Whang IS, Olsson CA, McMahon DJ, Cooner WH. "The Use of Prostate Specific Antigen Density to Enhance the Predictive Value of Intermediate Levels of Serum Prostate Antigen." *Journal of Urology* 147(3, Part 2):187, 1992.

Sections in the Patient Report

Section Name	Description
Exam Specific Information	Indications and additional information for specific exam types
2D Measurements	Measurement labels and results for 2D-mode measurements
M Measurements	Measurement labels and results for M-mode measurements
D Measurements	Measurement labels and results for Doppler measurements
Assessment	Descriptions for exam-related structures. You can select system defined or user-defined descriptions.
Prostate Assessment	Includes system-defined descriptions and observations for prostate structures.
Images	Images included from the thumbnail panel
Comments	Any combination of predefined comments and user-entered comments

Musculoskeletal Measurements and Calculations

The musculoskeletal exam type has only general measurements and no measurement labels.

Sections in the Patient Report

Section Name	Description
Exam Specific Information	Indications and additional information for specific exam types including systole and diastole for both left and right sides and the blood pressure for each side
2D Measurements	Measurement labels and results for 2D-mode measurements
2D Mode Side by Side	Measurement labels and results for 2D-mode measurements for both right and left sides
M Measurements	Measurement labels and results for M-mode measurements
D Measurements	Measurement labels and results for Doppler measurements
D Mode Side by Side	Measurement labels and results for Doppler measurements for both right and left sides
Assessment	Descriptions for exam-related structures. You can select system defined or user-defined descriptions.
Images	Images included from the thumbnail panel
Comments	Any combination of predefined comments and user-entered comments

Intraoperative Measurements and Calculations

The intraoperative exam type has only general measurements and no measurement labels.

Sections in the Patient Report

Section Name	Description
Exam Specific Information	Indications and additional information for specific exam types
2D Measurements	Measurement labels and results for 2D-mode measurements
2D Mode Side by Side	Measurement labels and results for 2D-mode measurements for both right and left sides
M Measurements	Measurement labels and results for M-mode measurements
D Measurements	Measurement labels and results for Doppler measurements
D Mode Side by Side	Measurement labels and results for Doppler measurements for both right and left sides
Assessment	Descriptions for exam-related structures. You can select system defined or user-defined descriptions.
Images	Images included from the thumbnail panel
Comments	Any combination of predefined comments and user-entered comments

Small Parts Measurements and Calculations

2D-Mode Measurement Labels

Measurement Label	Description	Measurement Method
Urethra	Urethra	Distance
Brachial A	Brachial Artery	Distance
Testicle Length	Testicle Length	Distance
Testicle AP	Testicle AP	Distance
Testicle Width	Testicle Width	Distance
Epididymis Length	Epididymis Length	Distance
Epididymis AP	Epididymis AP	Distance
Epididymis Width	Epididymis Width	Distance
Mass 1	Testicle Mass, numbered 1 through 5	3 Dist
Testicular A	Testicular Artery	Distance
Pampiniform Plexus	Pampiniform Plexus	Distance
Scrotal Wall	Scrotal Wall	Distance
Corp Cav	Corpus Cavernosom	Distance
Corp Spong	Corpus Spongiosum	Distance
Cav A	Cavernosal Artery	Distance
Pre-Inj Cav A	Pre-Injection Cavernosal Artery	Distance
Post-Inj Cav A	Post-Injection Cavernosal Artery	Distance
Dorsal A	Dorsal Artery	Distance
Pre-Inj Dorsal A	Pre-Injection Dorsal Artery	Distance
Post-Inj Dorsal A	Post-Injection Dorsal Artery	Distance

Measurement Label	Description	Measurement Method
Urethral A	Urethral Artery	Distance
Bulbar A	Bulbar Artery	Distance
Sup Dorsal V	Superficial Dorsal Vein	Distance
Dp Penile V	Deep Penile Vein	Distance
Penile V	Penile Vein	Distance
Isthmus AP	Anteroposterior diameter of thyroid isthmus	Distance
Thyroid Lobe TRV	Transverse diameter of thyroid lobe	Distance
Thyroid Lobe AP	Anteroposterior diameter of thyroid lobe	Distance
Thyroid Lobe SAG	Sagittal diameter of thyroid lobe	Distance
NOD 1	Thyroid nodule, numbered 1 through 15	3 Dist
Parathyroid	Parathyroid Length	Distance
Mass 1	Thyroid Mass, numbered 1 through 5	Distance
Isthmus Mass 1	Isthmus Mass, numbered 1 through 5	Distance
Neck Mass 1	Neck Mass, numbered 1 through 5	Distance
Superior Thyroid A	Superior Thyroid Artery	Distance
Inferior Thyroid A	Inferior Thyroid Artery	Distance

2D-Mode Calculation Labels

Calculation Label	Description	Required Measurements	Units
Thyroid Lobe Volume	Thyroid Lobe Volume $0.479 * \text{Thyroid Lobe SAG} * \text{Thyroid Lobe AP} * \text{Thyroid Lobe TRV}$	Thyroid Lobe SAG Thyroid Lobe TRV Thyroid Lobe AP	mL
Testicle Volume	Testicle Volume $3.1415926535 / 6 * \text{Testicle Length} * \text{Testicle Width} * \text{Testicle AP}$	Testicle Length Testicle Width Testicle AP	mL
Epididymis Volume	Epididymis Volume	Epididymis Length Epididymis Width Epididymis AP	mL

Thyroid Lobe Volume Reference

Brunn J, Block U, Ruf G, Bos I, Kunze WP, Scriba PC. "Volumetrie der Schilddrüsenlappen mittels Realtime-sonographie" *Deutsche Medizinische Wochenschrift* 106(41):1338–1340, 1981.

Doppler Measurement Labels

Measurement Label	Description	Measurement Method
Testicular V	Testicular Vein	Velocity
Brachial A	Brachial Artery	PI Auto
Testicular A	Testicular Artery	Velocity
Epididymal V	Epididymal Vein	Velocity
Intratesticular V	Intratesticular Vein	Velocity
Epididymal A	Epididymal Artery	RI
Intratesticular A	Intratesticular Artery	Velocity
Cav A	Cavernosal Artery	Velocity
Pre-Inj Cav A	Pre-Injection Cavernosal Artery	Velocity
Post-Inj Cav A	Post-Injection Cavernosal Artery	Velocity
Dorsal A	Dorsal Artery	Velocity
Pre-Inj Dorsal A	Pre-Injection Dorsal Artery	Velocity
Post-Inj Dorsal A	Post-Injection Dorsal Artery	Velocity
Urethral A	Urethral Artery	Velocity
Bulbar A	Bulbar Artery	Velocity
Sup Dorsal V	Superficial Dorsal Vein	Velocity
Dp Penile V	Deep Penile Vein	Velocity
Penile V	Penile Vein	Velocity
Cav A Vel	Cavernosal Artery Velocity	Velocity
Superior Thyroid A	Superior Thyroid Artery	PI Auto
Inferior Thyroid A	Inferior Thyroid Artery	PI Auto
Superior Thyroid V	Superior Thyroid Vein	Velocity
Middle Thyroid V	Middle Thyroid Vein	Velocity
Inferior Thyroid V	Inferior Thyroid Vein	Velocity

M-Mode Calculation Labels

Calculation Label	Description	Required Measurements	Units
LC / LBI	Ratio of the Left Cavernosal Artery to the Left Brachial Artery Cav A Vel / Brach A Vel	Cav A Vel Brach A Vel	n/a
RC / LBI	Ratio of the Right Cavernosal Artery to the Left Brachial Artery Cav A Vel / Brach A Vel	Cav A Vel Brach A Vel	n/a
LC / RBI	Ratio of the Left Cavernosal Artery to the Right Brachial Artery Cav A Vel / Brach A Vel	Cav A Vel Brach A Vel	n/a
RC / RBI	Ratio of the Right Cavernosal Artery to the Right Brachial Artery Cav A Vel / Brach A Vel	Cav A Vel Brach A Vel	n/a

Sections in the Patient Report

Section Name	Description
Exam Specific Information	Indications and additional information for specific exam types
2D Measurements	Measurement labels and results for 2D-mode measurements
2D Mode Side by Side	Measurement labels and results for 2D-mode measurements for both right and left sides
Assessment	Descriptions for exam-related structures. You can select system defined or user-defined descriptions.
Images	Images included from the thumbnail panel
Comments	Any combination of predefined comments and user-entered comments
M Measurements	Measurement labels and results for M-mode measurements
D Measurements	Measurement labels and results for Doppler measurements
D Mode Side by Side	Measurement labels and results for Doppler measurements for both right and left sides

Appendix F Obstetrical Bibliography

Obstetrical References.....	3
Mean Gestational Sac Diameter.....	3
Gestational Sac	3
Crown Rump Length.....	3
Biparietal Diameter	4
Occipital Frontal Diameter	4
Transverse Abdominal Diameter	4
Anterior-Posterior Abdominal Diameter	4
Fractional Shortening	4
Fetal Weight Estimation.....	5
Head Circumference.....	5
Abdominal Circumference	5
Anterior-Posterior Trunk Diameter x Transverse Trunk Diameter.....	5
Fetal Trunk Area.....	5
Femur Length	6
Humerus Length	6
Ulna Length	6
Binocular Distance.....	6
Tibia Length.....	6
Foot Length	6
Growth Analysis Ratios and Indexes.....	7
Uterine Artery: Pulsatility Index, Resistive Index	7
Cephalic Index.....	7
Fetal Weight Estimation.....	7
Femur Length/Abdominal Circumference Ratio.....	7
Femur Length/Biparietal Diameter Ratio.....	7
Head Circumference/Abdominal Circumference Ratio	7
Gestational Age by Ultrasound and Estimated Date of Delivery.....	8
Abdominal Circumference	8
Biparietal Diameter	8
Estimated Date of Delivery by Last Menstrual Period.....	8
Clinical Gestational Age.....	8
Crown Rump Length.....	8
Femur Length	9
Gestational Sac	9
Head Circumference.....	9
Humerus Length	9
Occipital Frontal Diameter	9
Cerebellum, Transverse Cerebellar Diameter	9
Composite Gestational Age Estimation by Ultrasound.....	10
Simple Average - Ultrasound Gestational Age (average)	10
Regression Equations, Hadlock.....	10
Estimated Date of Delivery by Ultrasound	10

Parameters for Growth Analysis Graphs.....	10
Mean Gestational Sac Diameter.....	10
Gestational Sac	10
Crown Rump Length.....	10
Biparietal Diameter	11
Biparietal Diameter (Outer to Inner).....	12
Biparietal Diameter (Outer to Outer).....	12
Occipital Frontal Diameter	12
Head Circumference.....	12
Head Circumference (Plotted)	13
Abdominal Circumference	13
Abdominal Circumference (Plotted).....	13
Anterior-Posterior Trunk Diameter x Transverse Trunk Diameter.....	13
Facial Angle.....	14
Fetal Trunk Area	14
Femur Length	14
Humerus Length	15
Ulna Length	15
Tibia Length.....	15
Foot Length	15
Fibula.....	15
Thoracic Circumference.....	16
Clavicle Length	16
Kidney Length.....	16
Kidney Anterior Posterior.....	16
Kidney Width	16
Transcerebellar Diameter	16
Transverse Abdominal Diameter	16
Cisterna Magna (Posterior Fossa).....	17
Anterior Cerebral Ventricle Diameter	17
Posterior Cerebral Ventricle Diameter	17
Umbilical Artery (Resistive Index).....	17
Umbilical Artery (Pulsatility Index)	17
Middle Cerebral Artery (Resistive Index)	17
Middle Cerebral Artery (Pulsatility Index).....	17
Mean Abdominal Diameter	17
Posterior Horn Lateral Ventricular Width	17
Cerebellum, Transverse Cerebellar Diameter	17
Estimated Fetal Weight for Growth Analysis Graphs.....	18
Radius	18
Ratios and Indices for Growth Analysis Graphs.....	18
Amniotic Fluid Index	18
Lateral Ventricular Width/Hemispheric Width Ratio	18
Cerebellum/Abdominal Circumference Ratio.....	18
Other Calculations.....	18
Corrected Biparietal Diameter	18

Obstetrical References

Mean Gestational Sac Diameter

Hellman LM, Kobayashi M, Fillisti L, Lavenhar M, Cromb E. "Growth and Development of the Human Fetus Prior to the Twentieth Week of Gestation." *American Journal of Obstetrics and Gynecology* 103(6):789, 1969.

Rempen A. "Biometrie in der Frühgravidität (I. Trimenon) (Biometry in Early Pregnancy (1st Trimester))." *Der Frauenarzt* 32:425, 1991.

Gestational Sac

Tokyo:

Masuda H, Shinozuka N, Okai T, Mizuno M. "Diagnosis of the Week of Pregnancy and Prognosis." *Perinatal Care* 8:719-726.

Crown Rump Length

Hadlock FP, Shah YP, Kanon DJ, Lindsey JV. "Fetal Crown-Rump Length: Reevaluation of Relation to Menstrual Age (5-18 weeks) with High-Resolution Real-Time US." *Radiology* 182(2):501, 1992.

Hansmann M, Hackelöer B-J, Staudach A. *Ultrasound Diagnosis in Obstetrics and Gynecology*. New York: Springer-Verlag, 1985, p. 439.

Lasser DM, Peisner DB, Vollebergh J, Timor-Tritsch I. "First-trimester fetal biometry using transvaginal sonography." *Ultrasound in Obstetrics and Gynecology* 3:104, 1993.

Robinson HP and Fleming JEE. "A critical evaluation of sonar 'crown-rump length' measurements." *British Journal of Obstetrics and Gynaecology* 82:702, 1975.

ASUM:

Westerway SC. "Ultrasonic Fetal Measurements: New Australian Standards for the New Millennium." *Aust NZ J Obstet Gynaecol* 40:3:297-302, 2000.

JSUM:

Japan Society of Ultrasonics in Medicine. "Standardization of Fetometry and Official Announcement of Diagnostic Guidelines." *J. Med. Ultrasonics* 30:415-440, 2003.

Osaka:

Aoki M, Yamada M. "Examining Fetal Growth." *Obstetrics and Gynecology* 47:547-556, 1983.

Tokyo:

Masuda H, Shinozuka N, Okai T, Mizuno M. "Diagnosis of the Week of Pregnancy and Prognosis." *Perinatal Care* 8:719-726.

Biparietal Diameter

Hadlock FP, Deter RL, Harrist RB, Park SK. "Estimating Fetal Age: Computer-Assisted Analysis of Multiple Fetal Growth Parameters." *Radiology* 152:497, 1984.

Lasser DM, Peisner DB, Vollebergh J, Timor-Tritsch I. "First-trimester fetal biometry using transvaginal sonography." *Ultrasound in Obstetrics and Gynecology* 3:104, 1993.

Merz E. *Ultrasound in Gynecology and Obstetrics*. Stuttgart and New York: Thieme Medical Publishers, Inc., 1991, p. 326.

Rempen A. "Biometrie in der Frühgravidität (I. Trimenon) (Biometry in Early Pregnancy (1st Trimester))." *Der Frauenarzt* 32:425, 1991.

ASUM:

Westerway SC. "Ultrasonic Fetal Measurements: New Australian Standards for the New Millennium." *Aust NZ J Obstet Gynaecol* 40:3:297-302, 2000.

JSUM:

Japan Society of Ultrasonics in Medicine. "Standardization of Fetometry and Official Announcement of Diagnostic Guidelines." *J. Med. Ultrasonics* 30:415-440, 2003.

Osaka:

Aoki M, Yamada M. "Examining Fetal Growth." *Obstetrics and Gynecology* 47:547-556, 1983.

Tokyo:

Masuda H, Shinozuka N, Okai T, Mizuno M. "Diagnosis of the Week of Pregnancy and Prognosis." *Perinatal Care* 8:719-726.

Occipital Frontal Diameter

ASUM

Westerway SC. "Ultrasonic Fetal Measurements: New Australian Standards for the New Millennium." *Aust NZ J Obstet Gynaecol* 40:3:297-302, 2000.

Transverse Abdominal Diameter

Merz E, Wellek S. "Normal Fetal Growth Profile – A Uniform Model for Calculating Normal Current Head and Abdomen Parameters and Long Limb Bones." *Ultraschall in der Medizin* 17(4):153-162, 1996.

Anterior-Posterior Abdominal Diameter

Merz E, Wellek S. "Normal Fetal Growth Profile – A Uniform Model for Calculating Normal Current Head and Abdomen Parameters and Long Limb Bones." *Ultraschall in der Medizin* 17(4):153-162, 1996.

Fractional Shortening

Cyr DR, Guntheroth WG, Mack LA. "Fetal Echocardiography." In: Berman MC (ed.) *Diagnostic Medical Sonography, Vol I: Obstetrics and Gynecology*, 249-271, 1991.

Fetal Weight Estimation

JSUM:

Japan Society of Ultrasonics in Medicine. "Standardization of Fetometry and Official Announcement of Diagnostic Guidelines." *J. Med. Ultrasonics* 30:415-440, 2003.

Osaka:

Aoki M, Yamada M. "Examining Fetal Growth." *Obstetrics and Gynecology* 47:547-556, 1983.

Tokyo:

Masuda H, Shinozuka N, Okai T, Mizuno M. "Diagnosis of the Week of Pregnancy and Prognosis." *Perinatal Care* 8:719-726.

Head Circumference

Hadlock FP, Deter RL, Harrist RB, Park SK. "Estimating Fetal Age: Computer-Assisted Analysis of Multiple Fetal Growth Parameters." *Radiology* 152:497, 1984.

Lasser DM, Peisner DB, Vollebergh J, Timor-Tritsch I. "First-trimester fetal biometry using transvaginal sonography." *Ultrasound in Obstetrics and Gynecology* 3:104, 1993.

Merz E. *Ultrasound in Gynecology and Obstetrics*. Stuttgart and New York: Thieme Medical Publishers, Inc., 1991, p. 326.

Abdominal Circumference

Hadlock FP, Deter RL, Harrist RB, Park SK. "Estimating Fetal Age: Computer-Assisted Analysis of Multiple Fetal Growth Parameters." *Radiology* 152:497, 1984.

Lasser DM, Peisner DB, Vollebergh J, Timor-Tritsch I. "First-trimester fetal biometry using transvaginal sonography." *Ultrasound in Obstetrics and Gynecology* 3:104, 1993.

Merz E. *Ultrasound in Gynecology and Obstetrics*. Stuttgart and New York: Thieme Medical Publishers, Inc., 1991, p. 326.

JSUM:

Japan Society of Ultrasonics in Medicine. "Standardization of Fetometry and Official Announcement of Diagnostic Guidelines." *J. Med. Ultrasonics* 30:415-440, 2003.

Anterior-Posterior Trunk Diameter x Transverse Trunk Diameter

JSUM:

Japan Society of Ultrasonics in Medicine. "Standardization of Fetometry and Official Announcement of Diagnostic Guidelines." *J. Med. Ultrasonics* 30:415-440, 2003.

Tokyo:

Masuda H, Shinozuka N, Okai T, Mizuno M. "Diagnosis of the Week of Pregnancy and Prognosis." *Perinatal Care* 8:719-726.

Fetal Trunk Area

Osaka:

Aoki M, Yamada M. "Examining Fetal Growth." *Obstetrics and Gynecology* 47:547-556, 1983.

Femur Length

Hadlock FP, Deter RL, Harrist RB, Park SK. "Estimating Fetal Age: Computer-Assisted Analysis of Multiple Fetal Growth Parameters." *Radiology* 152:497, 1984.

Jeanty P, Rodesch F, Delbeke D, Dumont JE. "Estimation of Gestational Age from Measurements of Fetal Long Bones." *Journal of Ultrasound in Medicine* 3:75, 1984.

Merz E. *Ultrasound in Gynecology and Obstetrics*. Stuttgart and New York: Thieme Medical Publishers, Inc., 1991, p. 326.

JSUM:

Japan Society of Ultrasonics in Medicine. "Standardization of Fetometry and Official Announcement of Diagnostic Guidelines." *J. Med. Ultrasonics* 30:415-440, 2003.

Osaka:

Aoki M, Yamada M. "Examining Fetal Growth." *Obstetrics and Gynecology* 47:547-556, 1983.

Tokyo:

Masuda H, Shinozuka N, Okai T, Mizuno M. "Diagnosis of the Week of Pregnancy and Prognosis." *Perinatal Care* 8:719-726.

Humerus Length

Jeanty P, Rodesch F, Delbeke D, Dumont JE. "Estimation of Gestational Age from Measurements of Fetal Long Bones." *Journal of Ultrasound in Medicine* 3:75, 1984.

Osaka:

Aoki M, Yamada M. "Examining Fetal Growth." *Obstetrics and Gynecology* 47:547-556, 1983.

Ulna Length

Jeanty P, Rodesch F, Delbeke D, Dumont JE. "Estimation of Gestational Age from Measurements of Fetal Long Bones." *Journal of Ultrasound in Medicine* 3:75, 1984.

Binocular Distance

Jeanty P, Cantraine F, Cousaert E, Romero R, Hobbins JC. "The Binocular Distance: A New Way to Estimate Fetal Age." *Journal of Ultrasound in Medicine* 3:241, 1984.

Tongsong T, Wanapirak C, Jesadapornchai S, Tathayathikom E. "Fetal binocular distance as a predictor of menstrual age." *International Journal of Gynecology and Obstetrics* 38:87, 1992.

Tibia Length

Jeanty P, Rodesch F, Delbeke D, Dumont JE. "Estimation of Gestational Age from Measurements of Fetal Long Bones." *Journal of Ultrasound in Medicine* 3:75, 1984.

Foot Length

Mercer BM, Sklar S, Shariatmadar A, Gillieson MS, Dalton ME. "Fetal foot length as a predictor of gestational age." *American Journal of Obstetrics and Gynecology* 156(2):350, 1987.

Growth Analysis Ratios and Indexes

Uterine Artery: Pulsatility Index, Resistive Index

Merz E. *Ultrasound in Obstetrics and Gynecology, Volume 1: Obstetrics*. Stuttgart and New York: Thieme Medical Publishers, Inc., 2005, 2nd Edition (ISBN 9783131318824); p. 469-480, 61.

Schaffer H. "IPDS 1998 11th Congress of the International Perinatal Doppler Society." *Journal für Fertilität und Reproduktion*. July 2-4, 1998.

Cephalic Index

Chitty, LS, Altman, DG. "Charts of Fetal Size: 2. Head Measurements." *British Journal of Obstetrics & Gynaecology* 101:35-43, 1994.

Hadlock FP, Deter RL, Carpenter RJ, Park SK. "Estimating Fetal Age: Effect of Head Shape on BPD." *American Journal of Roentgenology* 137:83, 1981.

Fetal Weight Estimation

Hadlock FP, Harrist RB, Sharman RS, Deter RL, Park SK. "Estimation of fetal weight with the use of head, body, and femur measurements—A prospective study." *American Journal of Obstetrics and Gynecology* 151:333, 1985.

Hansmann M, Hackelöer B-J, Staudach A. *Ultrasound Diagnosis in Obstetrics and Gynecology*. New York: Springer-Verlag, 1986.

Merz E. *Ultrasound in Gynecology and Obstetrics*. Stuttgart and New York: Thieme Medical Publishers, Inc., 1991, p. 157.

Shepard MJ, Richards VA, Berkowitz RL, Warsof SL, Hobbins JC. "An evaluation of two equations for predicting fetal weight by ultrasound." *American Journal of Obstetrics and Gynecology*. 142(1):47, 1982.

Femur Length/Abdominal Circumference Ratio

Hadlock FP, Deter RL, Harrist RB, Roecker E, Park SK. "A Date-Independent Predictor of Intrauterine Growth Retardation: Femur Length/Abdominal Circumference Ratio." *American Journal of Roentgenology* 141: 979, 1983.

Femur Length/Biparietal Diameter Ratio

Hohler CW, Quetel TA. "Comparison of ultrasound femur length and biparietal diameter in late pregnancy." *American Journal of Obstetrics and Gynecology* 141(7):759, 1981.

Head Circumference/Abdominal Circumference Ratio

Campbell S, Thoms A. "Ultrasound measurement of the fetal head to abdomen circumference ratio in the assessment of growth retardation." *British Journal of Obstetrics and Gynaecology*, 84:165, 1977.

Snijders RJ, Nicolaides KH. "Fetal biometry at 14-40 weeks' gestation." *Ultrasound in Obstetrics and Gynecology* 4:34-48, 1994.

Gestational Age by Ultrasound and Estimated Date of Delivery

Abdominal Circumference

ASUM 2003:

Westerway SC. "Ultrasonic Fetal Measurements for an Australian Population." 2003.
Adapted from Westerway SC, Davison A, Cowell S. "Ultrasonic Fetal Measurements: New Australian Standards for the New Millennium." *Australian and New Zealand Journal of Obstetrics and Gynaecology* 40(3):297–302, 2000.

CFEF:

Créquat J, Duyme M, Brodaty G. "Biométrie 2000. Tables de croissance fœtale par le Collège Français d'Echographie Fœtale (CFEF) et l'Inserm U155." ("Biometry 2000. Fetal growth charts by the French College of fetal ultrasonography and the Inserm U155.") *Gynécologie Obstétrique and Fertilité* 28(6):435–445, 2000.

Biparietal Diameter

Hansmann M, Hackelöer B-J, Staudach A. *Ultrasound Diagnosis in Obstetrics and Gynecology*. New York: Springer-Verlag, 1986, English Edition (ISBN 0-387-15384-9); p.440–441.

Jeanty P, Romero R. *Obstetrical Ultrasound*. New York: McGraw-Hill Book Company, 1984, (ISBN 0-070-32319-4); p. 58.

CFEF:

Créquat J, Duyme M, Brodaty G. "Biométrie 2000. Tables de croissance fœtale par le Collège Français d'Echographie Fœtale (CFEF) et l'Inserm U155." ("Biometry 2000. Fetal growth charts by the French College of fetal ultrasonography and the Inserm U155.") *Gynécologie Obstétrique and Fertilité* 28(6):435–445, 2000.

Chitty (Outer to Inner):

Altman DG, Chitty LS. "New Charts for Ultrasound Dating of Pregnancy." *Ultrasound in Obstetrics and Gynecology* 10(3):174–191, 1997.

Chitty (Outer to Outer):

Altman DG, Chitty LS. "New Charts for Ultrasound Dating of Pregnancy." *Ultrasound in Obstetrics and Gynecology* 10(3):174–191, 1997.

Verbarg (Outer to Outer):

Verbarg BO, Steegers EA, De Ridder M, Snijders RJ, Smith E, Hofman A, Moll HA, Jaddoe VW, Witteman JC. "New Charts for Ultrasound Dating of Pregnancy and Assessment of Fetal Growth: Longitudinal Data From a Population-based Cohort Study." *Ultrasound in Obstetrics and Gynecology* 31(4):388–396, 2008.

Estimated Date of Delivery by Last Menstrual Period

280 days added to date of last menstrual period.

Clinical Gestational Age

Source data for Clinical GA is determined by the user during the Patient ID operation calculated weeks and days since last menstrual period.

Crown Rump Length

Rempen A, Chaoui R, Kozłowski P, Häusler M, Terinde R, Wisser J. "Standards zur Ultraschalluntersuchung in der Frühschwangerschaft." *Der Frauenarzt* 42:327–331, 2001.

Verbarg BO, Steegers EA, De Ridder M, Snijders RJ, Smith E, Hofman A, Moll HA, Jaddoe VW, Witteman JC. "New Charts for Ultrasound Dating of Pregnancy and Assessment of Fetal Growth: Longitudinal Data From a Population-based Cohort Study." *Ultrasound in Obstetrics and Gynecology* 31(4):388–396, 2008.

INTERGROWTH-21st:

Papageorgiou AT, Kennedy SH, Salomon LJ, Ohuma EO, et al. "International Standards for Early Fetal Size and Pregnancy Dating Based on Ultrasound Measurement of Crown-rump Length in the First Trimester of Pregnancy." *Ultrasound in Obstetrics and Gynecology* 44(6):641–648, 2014.

Femur Length

Chitty LS, Altman DG. "Charts of fetal size: limb bones." *British Journal of Obstetrics & Gynaecology* 109:919-929, 2002.

Warda AH, Deter RL, Rossavik IK, Carpenter RJ, Hadlock FP. "Fetal Femur Length: A Critical Reevaluation of Relationship to Menstrual Age." *Obstetrics and Gynecology*, 66(1):69-75, 1985.

ASUM 2003:

Westerway SC. "Ultrasonic Fetal Measurements for an Australian Population." 2003.

Adapted from Westerway SC, Davison A, Cowell S. "Ultrasonic Fetal Measurements: New Australian Standards for the New Millennium." *Australian and New Zealand Journal of Obstetrics and Gynaecology* 40(3):297-302, 2000.

CFEF:

Créquat J, Duyme M, Brodaty G. "Biométrie 2000. Tables de croissance fœtale par le Collège Français d'Echographie Fœtale (CFEF) et l'Inserm U155." ("Biometry 2000. Fetal growth charts by the French College of fetal ultrasonography and the Inserm U155.") *Gynécologie Obstétrique and Fertilité* 28(6):435-445, 2000.

Gestational Sac

Rempen A. "Biometrie in der Frühgravidität (I. Trimenon) (Biometry in Early Pregnancy (1st Trimester))." *Der Frauenarzt* 32:425, 1991.

Head Circumference

Chitty, LS, Altman, DG. "Charts of Fetal Size: 2. Head Measurements." *British Journal of Obstetrics & Gynaecology* 101:35-43, 1994.

Hansmann M, Hackelöer B-J, Staudach A. *Ultrasound Diagnosis in Obstetrics and Gynecology*. New York: Springer-Verlag, 1986, English Edition (ISBN 0-387-15384-9); p. 440-441.

ASUM 2003:

Westerway SC. "Ultrasonic Fetal Measurements for an Australian Population." 2003.

Adapted from Westerway SC, Davison A, Cowell S. "Ultrasonic Fetal Measurements: New Australian Standards for the New Millennium." *Australian and New Zealand Journal of Obstetrics and Gynaecology* 40(3):297-302, 2000.

CFEF:

Créquat J, Duyme M, Brodaty G. "Biométrie 2000. Tables de croissance fœtale par le Collège Français d'Echographie Fœtale (CFEF) et l'Inserm U155." ("Biometry 2000. Fetal growth charts by the French College of fetal ultrasonography and the Inserm U155.") *Gynécologie Obstétrique and Fertilité* 28(6):435-445, 2000.

Humerus Length

ASUM 2003:

Westerway SC. "Ultrasonic Fetal Measurements for an Australian Population." 2003.

Adapted from Westerway SC, Davison A, Cowell S. "Ultrasonic Fetal Measurements: New Australian Standards for the New Millennium." *Australian and New Zealand Journal of Obstetrics and Gynaecology* 40(3):297-302, 2000.

Occipital Frontal Diameter

Chitty LS, Altman DG. "Charts of Fetal Size: 2. Head Measurements." *British Journal of Obstetrics and Gynaecology* 101(1):35-43, 1994.

Cerebellum, Transverse Cerebellar Diameter

Hill LM, Guzick D, Fries J, Hixson J, Rivello D. "The Transverse Cerebellar Diameter in Estimating Gestational Age in the Large for Gestational Age-Fetus." *Obstetrics and Gynecology* 75(6):981-985, 1990.

Composite Gestational Age Estimation by Ultrasound

Simple Average - Ultrasound Gestational Age (average)

The sum of all gestational ages determined by single parameter tables/formulas, divided by the number of parameters summed.

Regression Equations, Hadlock

Hadlock FP, Deter RL, Harrist RB, Park SK. "Estimating Fetal Age: Computer-Assisted Analysis of Multiple Fetal Growth Parameters." *Radiology* 152:497, 1984.

Estimated Date of Delivery by Ultrasound

280 days - today's ultrasound gestational age (in days), add to today's date = estimated date of delivery by ultrasound

Parameters for Growth Analysis Graphs

The references in this section are used to analyze fetal growth by plotting a measured parameter against clinical gestational age. The plotted points are displayed on a graph that indicates the 5%, 50%, and 95% for an expected parameter measurement determined by a function using gestational age.

Mean Gestational Sac Diameter

Hellman LM, Kobayashi M, Fillisti L, Lavenhar M, Cromb E. "Growth and Development of the Human Fetus Prior to the Twentieth Week of Gestation." *American Journal of Obstetrics and Gynecology* 103(6):789, 1969.

Rempen A. "Biometrie in der Frühgravidität (I. Trimenon) (Biometry in Early Pregnancy (1st Trimester))." *Der Frauenarzt* 32:425, 1991.

Gestational Sac

Hellman LM, Kobayashi M, Fillisti L, Lavenhar M, Cromb E. "Growth and Development of the Human Fetus Prior to the Twentieth Week of Gestation." *American Journal of Obstetrics and Gynecology* 103(6):789, 1969.

Rempen A. "Biometrie in der Frühgravidität (I. Trimenon) (Biometry in Early Pregnancy (1st Trimester))." *Der Frauenarzt* 32:425, 1991.

Tokyo:

Masuda H, Shinozuka N, Okai T, Mizuno M. "Diagnosis of the Week of Pregnancy and Prognosis." *Perinatal Care* 8:719-726.

Crown Rump Length

Hadlock FP, Shah YP, Kanon DJ, Lindsey JV. "Fetal Crown-Rump Length: Reevaluation of Relation to Menstrual Age (5-18 weeks) with High-Resolution Real-Time US." *Radiology* 182(2):501, 1992.

Hansmann M, Hackelöer B-J, Staudach A. *Ultrasound Diagnosis in Obstetrics and Gynecology*. New York: Springer-Verlag, 1985.

Pexsters A, Daemen A, Bottomley C, Van Schoubroeck D, De Catte L, et al. "New Crown-rump Length Curve Based on Over 3500 Pregnancies." *Ultrasound in Obstetrics and Gynecology* 35(6):650-655, 2010.

Rempen A. "Biometrie in der Frühgravidität (I. Trimenon) (Biometry in Early Pregnancy (1st Trimester))." *Der Frauenarzt* 32:425, 1991, p. 427.

Robinson HP and Fleming JEE. "A critical evaluation of sonar 'crown-rump length' measurements." *British Journal of Obstetrics and Gynaecology* 82:702, 1975.

ASUM:

Westerway SC. "Ultrasonic Fetal Measurements: New Australian Standards for the New Millennium." *Aust NZ J Obstet Gynaecol* 40:3:297-302, 2000.

INTERGROWTH-21st:

Papageorgiou AT, Kennedy SH, Salomon LJ, Ohuma EO, et al. "International Standards for Early Fetal Size and Pregnancy Dating Based on Ultrasound Measurement of Crown-rump Length in the First Trimester of Pregnancy." *Ultrasound in Obstetrics and Gynecology* 44(6):641–648, 2014.

JSUM:

Japan Society of Ultrasonics in Medicine. "Standardization of Fetometry and Official Announcement of Diagnostic Guidelines." *J. Med. Ultrasonics* 30:415-440, 2003.

Osaka:

Aoki M, Yamada M. "Examining Fetal Growth." *Obstetrics and Gynecology* 47:547-556, 1983.

Tokyo:

Masuda H, Shinozuka N, Okai T, Mizuno M. "Diagnosis of the Week of Pregnancy and Prognosis." *Perinatal Care* 8:719-726.

Biparietal Diameter

Hadlock FP, Deter RL, Harrist RB, Park SK. "Estimating Fetal Age: Computer-Assisted Analysis of Multiple Fetal Growth Parameters." *Radiology* 152:497, 1984.

Hansmann M, Hackelöer B-J, Staudach A. *Ultrasound Diagnosis in Obstetrics and Gynecology*. New York: Springer-Verlag, 1985, p. 176.

Lasser DM, Peisner DB, Vollebergh J, Timor-Tritsch I. "First-trimester fetal biometry using transvaginal sonography." *Ultrasound in Obstetrics and Gynecology* 3:104, 1993.

Merz E, Kim-Kern M-S, Pehl S. "Ultrasonic Mensuration of Fetal Limb Bones in the Second and Third Trimesters." *Journal of Clinical Ultrasound* 15:175, March/April 1987.

Rempen A. "Biometrie in der Frühgravidität (I. Trimenon) (Biometry in Early Pregnancy (1st Trimester))." *Der Frauenarzt* 32:425, 1991, p. 427.

Snijders RJ, Nicolaides KH. "Fetal biometry at 14-40 weeks' gestation." *Ultrasound in Obstetrics and Gynecology* 4:34-48, 1994.

Verburg BO, Steegers EA, De Ridder M, Snijders RJ, Smith E, Hofman A, et al. "New charts for ultrasound dating of pregnancy and assessment of fetal growth: longitudinal data from a population-based cohort study." *Ultrasound in Obstetrics and Gynecology* 31:388-396, 2008.

ASUM:

Westerway SC. "Ultrasonic Fetal Measurements: New Australian Standards for the New Millennium." *Aust NZ J Obstet Gynaecol* 40:3:297-302, 2000.

CFEF:

Créquat J, Duyme M, Brodaty G. "Biométrie 2000. Tables de croissance fœtale par le Collège Français d'Echographie Fœtale (CFEF) et l'Inserm U155." ("Biometry 2000. Fetal growth charts by the French College of fetal ultrasonography and the Inserm U155.") *Gynécologie Obstétrique and Fertilité* 28(6):435–445, 2000.

Grumbach (Twin):

Grumbach K, Coleman BG, Arger PH, Mintz MC, Gabbe SV, Mennuti MT. "Twin and Singleton Growth Patterns Compared Using US." *Radiology* 158(1):237–241, 1986.

INTERGROWTH-21st:

Papageorgiou AT, Ohuma EO, Altman DG, Todros T, Ismail LC, et al. "International Standards for Fetal Growth Based on Serial Ultrasound Measurements: The Fetal Growth Longitudinal Study of the INTERGROWTH-21st Project." *The Lancet* 384(9946):869–879, 2014.

JSUM:

Japan Society of Ultrasonics in Medicine. "Standardization of Fetometry and Official Announcement of Diagnostic Guidelines." *J. Med. Ultrasonics* 30:415-440, 2003.

Osaka:

Aoki M, Yamada M. "Examining Fetal Growth." *Obstetrics and Gynecology* 47:547-556, 1983.

Tokyo:

Masuda H, Shinozuka N, Okai T, Mizuno M. "Diagnosis of the Week of Pregnancy and Prognosis." *Perinatal Care* 8:719-726.

Biparietal Diameter (Outer to Inner)

Chitty, LS, Altman, DG. "Charts of Fetal Size: 2. Head Measurements." *British Journal of Obstetrics & Gynaecology* 101:35-43, 1994.

Biparietal Diameter (Outer to Outer)

Chitty, LS, Altman, DG. "Charts of Fetal Size: 2. Head Measurements." *British Journal of Obstetrics & Gynaecology* 101:35-43, 1994.

Occipital Frontal Diameter

Chitty, LS, Altman, DG. "Charts of Fetal Size: 2. Head Measurements." *British Journal of Obstetrics & Gynaecology* 101:35-43, 1994.

Hansmann M, Hackelöer B-J, Staudach A. *Ultrasound Diagnosis in Obstetrics and Gynecology*. New York: Springer-Verlag, 1985, p. 176.

Snijders RJ, Nicolaides KH. "Fetal biometry at 14-40 weeks' gestation." *Ultrasound in Obstetrics and Gynecology* 4:34-48, 1994.

ASUM:

Westerway SC. "Ultrasonic Fetal Measurements: New Australian Standards for the New Millennium." *Aust NZ J Obstet Gynaecol* 40:3:297-302, 2000.

INTERGROWTH-21st:

Papageorgiou AT, Ohuma EO, Altman DG, Todros T, Ismail LC, et al. "International Standards for Fetal Growth Based on Serial Ultrasound Measurements: The Fetal Growth Longitudinal Study of the INTERGROWTH-21st Project." *The Lancet* 384(9946):869-879, 2014.

Head Circumference

Deter RL, Harrist RB, Hadlock FP, Carpenter RJ. "Fetal Head and Abdominal Circumferences: II. A Critical Reevaluation of the Relationship to Menstrual Age." *Journal of Clinical Ultrasound* 10(8):365-372, 1982.

Hadlock FP, Deter RL, Harrist RB, Park SK. "Estimating Fetal Age: Computer-Assisted Analysis of Multiple Fetal Growth Parameters." *Radiology* 152:497, 1984.

Hansmann M, Hackelöer B-J, Staudach A. *Ultrasound Diagnosis in Obstetrics and Gynecology*. New York: Springer-Verlag, 1985, p. 176.

Merz E. *Ultrasound in Gynecology and Obstetrics*. Stuttgart and New York: Thieme Medical Publishers, Inc., 1991, p. 312.

Snijders RJ, Nicolaides KH. "Fetal biometry at 14-40 weeks' gestation." *Ultrasound in Obstetrics and Gynecology* 4:34-48, 1994.

Verburg BO, Steegers EA, De Ridder M, Snijders RJ, Smith E, Hofman A, et al. "New charts for ultrasound dating of pregnancy and assessment of fetal growth: longitudinal data from a population-based cohort study." *Ultrasound in Obstetrics and Gynecology* 31:388-396, 2008.

ASUM:

Westerway SC. "Ultrasonic Fetal Measurements: New Australian Standards for the New Millennium." *Aust NZ J Obstet Gynaecol* 40:3:297-302, 2000.

CFEF:

Créquat J, Duyme M, Brodaty G. "Biométrie 2000. Tables de Croissance Fœtale par le Collège Français d'Echographie Fœtale (CFEF) et l'Inserm U155." ("Biometry 2000. Fetal Growth Charts by the French College of Fetal Ultrasonography and the Inserm U155.") *Gynécologie Obstétrique and Fertilité* 28(6):435-445, 2000.

INTERGROWTH-21st:

Papageorgiou AT, Ohuma EO, Altman DG, Todros T, Ismail LC, et al. "International Standards for Fetal Growth Based on Serial Ultrasound Measurements: The Fetal Growth Longitudinal Study of the INTERGROWTH-21st Project." *The Lancet* 384(9946):869-879, 2014.

Head Circumference (Plotted)

Chitty, LS, Altman, DG. "Charts of Fetal Size: 2. Head Measurements." *British Journal of Obstetrics & Gynaecology* 101:35-43, 1994.

Abdominal Circumference

Deter RL, Harrist RB, Hadlock FP, Carpenter RJ. "Fetal Head and Abdominal Circumferences: II. A Critical Reevaluation of the Relationship to Menstrual Age." *Journal of Clinical Ultrasound* 10(8):365-372, 1982.

Hadlock FP, Deter RL, Harrist RB, Park SK. "Estimating Fetal Age: Computer-Assisted Analysis of Multiple Fetal Growth Parameters." *Radiology* 152:497, 1984.

Jeanty P, Cousaert E, Cantraine F. "Normal growth of the abdominal perimeter." *American Journal of Perinatology* 1:136, 1984, p. 133.

Merz E. *Ultrasound in Gynecology and Obstetrics*. Stuttgart and New York: Thieme Medical Publishers, Inc., 1991, p. 312.

Snijders RJ, Nicolaides KH. "Fetal biometry at 14-40 weeks' gestation." *Ultrasound in Obstetrics and Gynecology* 4:34-48, 1994.

Verburg BO, Steegers EA, De Ridder M, Snijders RJ, Smith E, Hofman A, et al. "New charts for ultrasound dating of pregnancy and assessment of fetal growth: longitudinal data from a population-based cohort study." *Ultrasound in Obstetrics and Gynecology* 31:388-396, 2008.

ASUM:

Westerway SC. "Ultrasonic Fetal Measurements: New Australian Standards for the New Millennium." *Aust NZ J Obstet Gynaecol* 40:3:297-302, 2000.

CFEF:

Créquat J, Duyme M, Brodaty G. "Biométrie 2000. Tables de croissance fœtale par le Collège Français d'Echographie Fœtale (CFEF) et l'Inserm U155." ("Biometry 2000. Fetal growth charts by the French College of fetal ultrasonography and the Inserm U155.") *Gynécologie Obstétrique and Fertilité* 28(6):435-445, 2000.

Grumbach (Twin):

Grumbach K, Coleman BG, Arger PH, Mintz MC, Gabbe SV, Mennuti MT. "Twin and Singleton Growth Patterns Compared Using US." *Radiology* 158(1):237-241, 1986.

INTERGROWTH-21st:

Papageorgiou AT, Ohuma EO, Altman DG, Todros T, Ismail LC, et al. "International Standards for Fetal Growth Based on Serial Ultrasound Measurements: The Fetal Growth Longitudinal Study of the INTERGROWTH-21st Project." *The Lancet* 384(9946):869-879, 2014.

JSUM:

Japan Society of Ultrasonics in Medicine. "Standardization of Fetometry and Official Announcement of Diagnostic Guidelines." *J. Med. Ultrasonics* 30:415-440, 2003.

Abdominal Circumference (Plotted)

Chitty, LS, Altman, DG. "Charts of Fetal Size: 2. Head Measurements." *British Journal of Obstetrics & Gynaecology* 101:35-43, 1994.

Anterior-Posterior Trunk Diameter x Transverse Trunk Diameter

Tokyo:

Masuda H, Shinozuka N, Okai T, Mizuno M. "Diagnosis of the Week of Pregnancy and Prognosis." *Perinatal Care* 8:719-726.

JSUM:

Japan Society of Ultrasonics in Medicine. "Standardization of Fetometry and Official Announcement of Diagnostic Guidelines." *J. Med. Ultrasonics* 30:415-440, 2003.

Facial Angle

Sonek J. "Fronto-Maxillary Facial Angle in Screening for Trisomy 21." *Fetal Medicine Foundation of the United States of America* 3(1):1-2, 2007.

Fetal Trunk Area

Osaka:

Aoki M, Yamada M. "Examining Fetal Growth." *Obstetrics and Gynecology* 47:547-556, 1983.

Femur Length

Chitty LS, Altman DG. "Charts of fetal size: limb bones." *British Journal of Obstetrics & Gynaecology* 109:919-929, August 2002.

Hadlock FP, Deter RL, Harrist RB, Park SK. "Estimating Fetal Age: Computer-Assisted Analysis of Multiple Fetal Growth Parameters." *Radiology* 152:497, 1984.

Hansmann M, Hackelöer B-J, Staudach A. *Ultrasound Diagnosis in Obstetrics and Gynecology*. New York: Springer-Verlag, 1985.

Jeanty P, Dramaix-Wilmet M, van Kerkem J, Petroons P, Schwes J. "Ultrasonic Evaluation of Fetal Limb Growth, Part II." *Radiology* 143:751, 1982.

Merz E, Kim-Kern M-S, Pehl S. "Ultrasonic Mensuration of Fetal Limb Bones in the Second and Third Trimesters." *Journal of Clinical Ultrasound* 15:175, March/April 1987.

Snijders RJ, Nicolaides KH. "Fetal biometry at 14-40 weeks' gestation." *Ultrasound in Obstetrics and Gynecology* 4:34-48, 1994.

Verburg BO, Steegers EA, De Ridder M, Snijders RJ, Smith E, Hofman A, et al. "New charts for ultrasound dating of pregnancy and assessment of fetal growth: longitudinal data from a population-based cohort study." *Ultrasound in Obstetrics and Gynecology* 31:388-396, 2008.

Warda AH, Deter RL, Rossavik IK, Carpenter RJ, Hadlock FP. "Fetal Femur Length: A Critical Reevaluation of Relationship to Menstrual Age." *Obstetrics and Gynecology* 66(1):69-75, 1985.

ASUM:

Australian Society of Ultrasound Medicine (ASUM) - Policy on Obstetric Exams, 1998.

CFEF:

Créquat J, Duyme M, Brodaty G. "Biométrie 2000. Tables de croissance fœtale par le Collège Français d'Echographie Fœtale (CFEF) et l'Inserm U155." ("Biometry 2000. Fetal growth charts by the French College of fetal ultrasonography and the Inserm U155.") *Gynécologie Obstétrique and Fertilité* 28(6):435-445, 2000.

INTERGROWTH-21st:

Papageorgiou AT, Ohuma EO, Altman DG, Todros T, Ismail LC, et al. "International Standards for Fetal Growth Based on Serial Ultrasound Measurements: The Fetal Growth Longitudinal Study of the INTERGROWTH-21st Project." *The Lancet* 384(9946):869-879, 2014.

JSUM:

Japan Society of Ultrasonics in Medicine. "Standardization of Fetometry and Official Announcement of Diagnostic Guidelines." *J. Med. Ultrasonics* 30:415-440, 2003.

Osaka:

Aoki M, Yamada M. "Examining Fetal Growth." *Obstetrics and Gynecology* 47:547-556, 1983.

Tokyo:

Masuda H, Shinozuka N, Okai T, Mizuno M. "Diagnosis of the Week of Pregnancy and Prognosis." *Perinatal Care* 8:719-726.

Humerus Length

Chitty LS, Altman DG. "Charts of fetal size: limb bones." *British Journal of Obstetrics & Gynaecology* 109:919-929, August 2002.

Hansmann M, Hackelöer B-J, Staudach A. *Ultrasound Diagnosis in Obstetrics and Gynecology*. New York: Springer-Verlag, 1985.

Jeanty P, Dramaix-Wilmet M, van Kerkem J, Petroons P, Schwes J. "Ultrasonic Evaluation of Fetal Limb Growth, Part II." *Radiology* 143:751, 1982.

Merz E, Kim-Kern M-S, Pehl S. "Ultrasonic Mensuration of Fetal Limb Bones in the Second and Third Trimesters." *Journal of Clinical Ultrasound* 15:175, March/April 1987.

ASUM:

Westerway SC. "Ultrasonic Fetal Measurements: New Australian Standards for the New Millennium." *Aust NZ J Obstet Gynaecol* 40:3:297-302, 2000.

Osaka:

Aoki M, Yamada M. "Examining Fetal Growth." *Obstetrics and Gynecology* 47:547-556, 1983.

Ulna Length

Chitty LS, Altman DG. "Charts of fetal size: limb bones." *British Journal of Obstetrics & Gynaecology* 109:919-929, August 2002.

Hansmann M, Hackelöer B-J, Staudach A. *Ultrasound Diagnosis in Obstetrics and Gynecology*. New York: Springer-Verlag, 1986.

Jeanty P, Dramaix-Wilmet M, van Kerkem J, Petroons P, Schwes J. "Ultrasonic Evaluation of Fetal Limb Growth, Part II." *Radiology* 143:751, 1982.

Merz E, Kim-Kern M-S, Pehl S. "Ultrasonic Mensuration of Fetal Limb Bones in the Second and Third Trimesters." *Journal of Clinical Ultrasound* 15:175, March/April 1987.

Tibia Length

Chitty LS, Altman DG. "Charts of fetal size: limb bones." *British Journal of Obstetrics & Gynaecology* 109:919-929, August 2002.

Hansmann M, Hackelöer B-J, Staudach A. *Ultrasound Diagnosis in Obstetrics and Gynecology*. New York: Springer-Verlag, 1985.

Jeanty P, Dramaix-Wilmet M, van Kerkem J, Petroons P, Schwes J. "Ultrasonic Evaluation of Fetal Limb Growth, Part II." *Radiology* 143:751, 1982.

Merz E, Kim-Kern M-S, Pehl S. "Ultrasonic Mensuration of Fetal Limb Bones in the Second and Third Trimesters." *Journal of Clinical Ultrasound* 15:175, March/April 1987.

Foot Length

Chitty LS, Altman DG. "Charts of fetal size: limb bones." *British Journal of Obstetrics & Gynaecology* 109:919-929, August 2002.

Mercer BM, Sklar S, Shariatmadar A, Gillieson MS, Dalton ME. "Fetal foot length as a predictor of gestational age." *American Journal of Obstetrics and Gynecology* 156(2):350, 1987.

Fibula

Chitty LS, Altman DG. "Charts of fetal size: limb bones." *British Journal of Obstetrics & Gynaecology* 109:919-929, 2002.

Thoracic Circumference

Chitkara U, Rosenberg J, Chervenak FA, Berkowitz GS, Levine R, Fagerstrom RM, Walker B, Berkowitz RL. "Prenatal sonographic assessment of the fetal thorax: Normal values." *American Journal of Obstetrics and Gynecology* 156:1069, 1987.

Clavicle Length

Yarkoni S, Schmidt W, Jeanty P, Reece EA, Hobbins JC. "Clavicular Measurement: A New Biometric Parameter for Fetal Evaluation." *Journal of Ultrasound in Medicine* 4:467, 1985.

Kidney Length

Bertagnoli L, Lalatta F, Gallicchio R, Fantuzzi M, Rusca M, Zorzoli A, Deter RL. "Quantitative Characterization of the Growth of the Fetal Kidney." *Journal of Clinical Ultrasound* 11:349, 1983.

Hansmann M, Hackelöer B-J, Staudach A. *Ultrasound Diagnosis in Obstetrics and Gynecology*. New York: Springer-Verlag, 1985.

CFEF:

Van Vuuren SH, Damen-Elias HA, Stigter RH, et al. "Size and Volume Charts of Fetal Kidney, Renal Pelvis and Adrenal Gland." *Ultrasound in Obstetrics and Gynecology* 40:659-664, 2012.

Kidney Anterior Posterior

Bertagnoli L, Lalatta F, Gallicchio R, Fantuzzi M, Rusca M, Zorzoli A, Deter RL. "Quantitative Characterization of the Growth of the Fetal Kidney." *Journal of Clinical Ultrasound* 11:349, 1983.

Hansmann M, Hackelöer B-J, Staudach A. *Ultrasound Diagnosis in Obstetrics and Gynecology*. New York: Springer-Verlag, 1985.

CFEF:

Van Vuuren SH, Damen-Elias HA, Stigter RH, et al. "Size and Volume Charts of Fetal Kidney, Renal Pelvis and Adrenal Gland." *Ultrasound in Obstetrics and Gynecology* 40(6):659-664, 2012.

Atoub F, Favre X. "Hauteur de la surrenale." Collège Français D'Échographie Fœtale. 2014. https://www.cfef.org/boite_a_outils/images/hauteursurrenale2014.pdf

Kidney Width

CFEF:

Van Vuuren SH, Damen-Elias HA, Stigter RH, et al. "Size and Volume Charts of Fetal Kidney, Renal Pelvis and Adrenal Gland." *Ultrasound in Obstetrics and Gynecology* 40(6):659-664, 2012.

Atoub F, Favre X. "Hauteur de la surrenale." Collège Français D'Échographie Fœtale. 2014. https://www.cfef.org/boite_a_outils/images/hauteursurrenale2014.pdf

Transcerebellar Diameter

Snijders RJ, Nicolaides KH. "Fetal biometry at 14-40 weeks' gestation." *Ultrasound in Obstetrics and Gynecology* 4:34-48, 1994.

Verbarg BO, Steegers EA, De Ridder M, Snijders RJ, Smith E, Hofman A, et al. "New charts for ultrasound dating of pregnancy and assessment of fetal growth: longitudinal data from a population-based cohort study." *Ultrasound in Obstetrics and Gynecology* 31:388-396, 2008.

Transverse Abdominal Diameter

CFEF:

Créquat J, Duyme M, Brodaty G. "Biométrie 2000. Tables de croissance fœtale par le Collège Français d'Echographie Fœtale (CFEF) et l'Inserm U155." ("Biometry 2000. Fetal growth charts by the French College of fetal ultrasonography and the Inserm U155.") *Gynécologie Obstétrique and Fertilité* 28(6):435-445, 2000.

Cisterna Magna (Posterior Fossa)

Snijders RJ, Nicolaides KH. "Fetal biometry at 14-40 weeks' gestation." *Ultrasound in Obstetrics and Gynecology* 4:34-48, 1994.

Anterior Cerebral Ventricle Diameter

Snijders RJ, Nicolaides KH. "Fetal biometry at 14-40 weeks' gestation." *Ultrasound in Obstetrics and Gynecology* 4:34-48, 1994.

Posterior Cerebral Ventricle Diameter

Snijders RJ, Nicolaides KH. "Fetal biometry at 14-40 weeks' gestation." *Ultrasound in Obstetrics and Gynecology* 4:34-48, 1994.

Umbilical Artery (Resistive Index)

JSUM:

Japan Society of Ultrasonics in Medicine. "Standardization of Fetometry and Official Announcement of Diagnostic Guidelines." *J. Med. Ultrasonics* 30:415-440, 2003.

Umbilical Artery (Pulsatility Index)

JSUM:

Japan Society of Ultrasonics in Medicine. "Standardization of Fetometry and Official Announcement of Diagnostic Guidelines." *J. Med. Ultrasonics* 28:844-872, 2001.

Middle Cerebral Artery (Resistive Index)

JSUM:

Japan Society of Ultrasonics in Medicine. "Standardization of Fetometry and Official Announcement of Diagnostic Guidelines." *J. Med. Ultrasonics* 30:415-440, 2003.

Middle Cerebral Artery (Pulsatility Index)

JSUM:

Japan Society of Ultrasonics in Medicine. "Standardization of Fetometry and Official Announcement of Diagnostic Guidelines." *J. Med. Ultrasonics* 30:415-440, 2003.

Mean Abdominal Diameter

Kurmanavicius J, Huch A, et al. "Standard – Tabellen und Kurven für Ultraschall- / Dopplersonografie." UniversitätsSpital Zürich, Dept. Frauenheilkunde Klinik und Poliklinik für Geburtshilfe; 1997.

Posterior Horn Lateral Ventricular Width

Snijders RJ, Nicolaides KH. "Fetal biometry at 14-40 weeks' gestation." *Ultrasound in Obstetrics and Gynecology* 4:34-48, 1994.

Cerebellum, Transverse Cerebellar Diameter

Snijders RJ, Nicolaides KH. "Fetal Biometry at 14–40 Weeks' Gestation." *Ultrasound in Obstetrics and Gynecology* 4(1):34–48, 1994.

Verburg BO, Steegers EA, De Ridder M, Snijders RJ, Smith E, Hofman A, Moll HA, Jaddoe VW, Witteman JC. "New Charts for Ultrasound Dating of Pregnancy and Assessment of Fetal Growth: Longitudinal Data From a Population-based Cohort Study." *Ultrasound in Obstetrics and Gynecology* 31(4):388–396, 2008.

CFEF:

Chavez MR, Ananth CV, Smulian JC, Lashley S, Kontopoulos EV, et al. "Fetal transcerebellar diameter nomogram in singleton gestations with special emphasis in the third trimester: a comparison with previously published nomograms." *American Journal of Obstetrics and Gynecology* 189:1021–1025, 2003.

Estimated Fetal Weight for Growth Analysis Graphs

- Brenner WE, Edelman DA, Hendricks CH. "A Standard of Fetal Growth for the United States of America." *American Journal of Obstetrics and Gynecology* 126(5):555–564, 1976.
- Doubilet PT, Benson CB, Nadel AS, Ringer SA. "Improved Birth Weight Table for Neonates Developed from Gestations Dated by Early Ultrasonography." *Journal of Ultrasound in Medicine* 16(4):241–249, 1997.
- Hadlock FP, Harrist RB, Martinez-Poyer J. "In Utero Analysis of Fetal Growth: A Sonographic Weight Standard." *Radiology* 181:129, 1991.
- Hansmann M, Hackelöer B-J, Staudach A. *Ultrasound Diagnosis in Obstetrics and Gynecology*. New York: Springer-Verlag, 1986.
- Jeanty P, Cantraine F, Romero R, Coussaert E, Hobbins JC. "A Longitudinal Study of Fetal Weight Growth." *Journal of Ultrasound in Medicine* 3:321, 1984.
- Williams RL, Creasy RK, Cunningham GC, Hawes WE, Norris FD, Tashiro M. "Fetal Growth and Perinatal Viability in California." *Obstetrics and Gynecology* 59:624:632, 1982.
- Yarkoni S, Reece EA, Holford T, O'Connor TZ, Hobbins JC. "Estimated Fetal Weight in the Evaluation of Growth in Twin Gestations: A Prospective Longitudinal Study." *Obstetrics and Gynecology* 69:636, 1987.

JSUM:

- Japan Society of Ultrasonics in Medicine. "Standardization of Fetometry and Official Announcement of Diagnostic Guidelines." *J. Med. Ultrasonics* 30:415-440, 2003.

Osaka:

- Aoki M, Yamada M. "Examining Fetal Growth." *Obstetrics and Gynecology* 47:547-556, 1983.

Tokyo:

- Masuda H, Shinozuka N, Okai T, Mizuno M. "Diagnosis of the Week of Pregnancy and Prognosis." *Perinatal Care* 8:719-726.

Radius

- Chitty LS, Altman DG. "Charts of fetal size: limb bones." *British Journal of Obstetrics & Gynaecology* 109:919-929, August 2002.

Ratios and Indices for Growth Analysis Graphs

Amniotic Fluid Index

- Moore TR, Cayle JE. "The amniotic fluid index in normal human pregnancy." *American Journal of Obstetrics and Gynecology* 162:1168, 1990.

Lateral Ventricular Width/Hemispheric Width Ratio

- Johnson ML, Dunne MG, Mack LA, Rashbaum CL. "Evaluation of Fetal Intracranial Anatomy by Static and Real-Time Ultrasound." *Journal of Clinical Ultrasound* 8:311, 1980.

Cerebellum/Abdominal Circumference Ratio

- Meyer WJ, Gauthier DW, Goldenberg B, Santolaya J, Sipos J, Cattledge F. "The Fetal Transverse Cerebellar Diameter/Abdominal Circumference Ratio: A Gestational Age-Independent Method of Assessing Fetal Size." *Journal of Ultrasound in Medicine* 12:379, 1993.
- Snijders RJ, Nicolaides KH. "Fetal biometry at 14-40 weeks' gestation." *Ultrasound in Obstetrics and Gynecology* 4:34-48, 1994.

Other Calculations

Corrected Biparietal Diameter

- Doubilet PM, Greenes RA. "Improved Prediction of Gestational Age from Fetal Head Measurements." *American Journal of Roentgenology*. 142:797, 1984.

Appendix G Cardiac Bibliography

Body Surface Area.....	3
2D-Mode Measurements and Calculations.....	3
Cardiac Index	3
Cardiac Output	3
Left Ventricular End-Diastolic Volume and Left Ventricular End-Systolic Volume	3
Teichholz Formula	3
Simpson Single Plane Formula	3
Simpson Bi-Plane.....	3
Ejection Fraction.....	3
Fractional Shortening	3
Stroke Index	3
Stroke Volume	4
Volume	4
M-Mode Measurements and Calculations	4
Cardiac Index	4
Cardiac Output	4
Left Ventricular End-Diastolic Volume and Left Ventricular End-Systolic Volume	4
Teichholz Formula.....	4
Ejection Fraction.....	4
Fractional Shortening	4
Aorta/Left Atrium Ratio	4
Heart Rate	5
Stroke Index	5
Stroke Volume.....	5
Doppler Mode Measurements and Calculations	5
Doppler References.....	5
Acceleration Slope/Time.....	5
Deceleration Slope/Time	5
Pressure Half-Time.....	5
Mitral Valve Area Continuity Equation	5
Aortic Valve Area.....	5
Heart Rate	5
Cardiac Output	6
Mitral Valve, E point to A point Ratio	6
Mitral Valve, A point to E point Ratio	6
Mitral Valve, CA amplitude to CE amplitude Ratio.....	6

Isovolumic Relaxation Time	6
Left Ventricular Index of Myocardial Performance	6
Right Ventricular Index of Myocardial Performance	6
Right Ventricular Systolic Pressure	6
Tricuspid Regurgitant Jet Plus Right Atrial Pressure	6
Maximum Pressure Gradient	7
Mean Pressure Gradient	7
Aortic Valve maximum Pressure Gradient	7
Aortic Valve mean Pressure Gradient	7
Mitral Valve maximum Pressure Gradient	7
Mitral Valve mean Pressure Gradient	7
Mitral Regurgitation maximum Pressure Gradient	7
Pulmonary Valve maximum Pressure Gradient	7
Pulmonary Valve mean Pressure Gradient	7
Tricuspid Regurgitation maximum Pressure Gradient	8
Tricuspid Regurgitation mean Pressure Gradient	8
Tricuspid Valve maximum Pressure Gradient	8
Tricuspid Valve mean Pressure Gradient	8
Left Ventricular Outflow Tract maximum Pressure Gradient	8
Left Ventricular Outflow Tract mean Pressure Gradient	8
Ventricular Septal Defect maximum Pressure Gradient	8
Additional Cardiac References	8
Left Atrial and Right Atrial Volumes	8
Left Ventricular Mass	9
Mitral Valve Area and Aortic Valve Area Using Continuity Equation with both the Peak Velocity and Velocity-time Integral	9
Cardiac Normal Range Values	10
Shortening / Wall Stress	10
Relative Wall Thickness	10
syngo VVI Clinical References	10

Body Surface Area

Du Bois D, Du Bois EF. "A Formula to Estimate the Approximate Surface Area if Height and Weight Be Known." *Archives of Internal Medicine* 17:863, 1916, as reprinted in *Nutrition* 5(5):303-311, 1989.

2D-Mode Measurements and Calculations

Cardiac Index

Berkow R, Editor. *The Merck Manual of Diagnosis and Therapy, 16th Edition*. Rahway, New Jersey: Merck and Co., 1992, p. 378.

Schiller NB, Shah PM, Crawford M, et al. "Recommendations for Quantitation of the Left Ventricle by Two-Dimensional Echocardiography." *Journal of the American Society of Echocardiography* 2(5):364, 1989.

Cardiac Output

Belenkie I, Nutter DO, Clark DW, McCraw DB, Raizner AE. "Assessment of Left Ventricular Dimensions and Function by Echocardiography." *American Journal of Cardiology* 31:755-762, June 1973.

Haïtes NE, McLennan FM, Mowat DHR, Rawles JM. "Assessment of Cardiac Output by Ultrasound Technique Alone." *British Heart Journal* 53:123-129, 1985.

Huntsman LL, Stewart DK, Barnes SR, Franklin SB, Colocousis JS, Hessel EA. "Noninvasive Doppler Determination of Cardiac Output in Man." *Circulation* 67(3):593-602, 1983.

Ihlen H, Amlie JP, Dale J, et. al. "Determination of cardiac output by Doppler echocardiography." *British Heart Journal* 51:54-60, 1984.

Left Ventricular End-Diastolic Volume and Left Ventricular End-Systolic Volume

Teichholz Formula

Teichholz LE, Kreulen T, Herman MV, Gorlin R. "Problems in Echocardiographic Volume Determinations: Echocardiographic-Angiographic Correlations in the Presence or Absence of Asynrgy." *American Journal of Cardiology* 37(1):7-11, 1976.

Simpson Single Plane Formula

Schiller NB, Shah PM, Crawford M, et al. "Recommendations for Quantitation of the Left Ventricle by Two-Dimensional Echocardiography." *Journal of the American Society of Echocardiography* 2(5):364, 1989.

Simpson Bi-Plane

Schiller NB, Shah PM, Crawford M, et al. "Recommendations for Quantitation of the Left Ventricle by Two-Dimensional Echocardiography." *Journal of the American Society of Echocardiography* 2(5):364, 1989.

Ejection Fraction

Pombo JF, Troy BL, Russell RO. "Left Ventricular Volumes and Ejection Fraction by Echocardiography." *Circulation* 43:480-490, 1971.

Fractional Shortening

Belenkie I, Nutter DO, Clark DW, McCraw DB, Raizner AE. "Assessment of Left Ventricular Dimensions and Function by Echocardiography." *American Journal of Cardiology* 31:755-762, June 1973.

Stroke Index

Görge G, Erbel R, Brennecke R, Rupprecht HJ, Todt M, Meyer J. "High-Resolution Two-dimensional Echocardiography Improves the Quantification of Left Ventricular Function." *Journal of the American Society of Echocardiography* 5(2):125-134, 1992.

Stroke Volume

Görge G, Erbel R, Brennecke R, Rupprecht HJ, Todt M, Meyer J. "High-Resolution Two-dimensional Echocardiography Improves the Quantification of Left Ventricular Function." *Journal of the American Society of Echocardiography* 5(2):125-134, 1992.

Roelandt J. *Practical Echocardiography* Ultrasound in Biomedicine, No. 1. Research Studies Press (a division of John Wiley & Sons Ltd), 1977, p. 270.

Volume

Schiller NB, Shah PM, Crawford M, et al. "Recommendations for Quantitation of the Left Ventricle by Two-Dimensional Echocardiography." *Journal of the American Society of Echocardiography* 2(5):364, 1989.

M-Mode Measurements and Calculations

Cardiac Index

Berkow, R, Editor. *The Merck Manual of Diagnosis and Therapy, 16th Edition*. Rahway, New Jersey: Merck and Co., 1992, p. 378.

Schiller NB, Shah PM, Crawford M, et al. "Recommendations for Quantitation of the Left Ventricle by Two-Dimensional Echocardiography." *Journal of the American Society of Echocardiography* 2(5):364, 1989.

Cardiac Output

Belenkie I, Nutter DO, Clark DW, McCraw DB, Raizner AE. "Assessment of Left Ventricular Dimensions and Function by Echocardiography." *American Journal of Cardiology* 31:755-762, June 1973.

Haïtes NE, McLennan FM, Mowat DHR, Rawles JM. "Assessment of Cardiac Output by Ultrasound Technique Alone." *British Heart Journal* 53:123-129, 1985.

Huntsman LL, Stewart DK, Barnes SR, Franklin SB, Colocousis JS, Hessel EA. "Noninvasive Doppler Determination of Cardiac Output in Man." *Circulation* 67(3):593-602, 1983.

Ihlen H, Amlie JP, Dale J, et. al. "Determination of cardiac output by Doppler echocardiography." *British Heart Journal* 51:54-60, 1984.

Left Ventricular End-Diastolic Volume and Left Ventricular End-Systolic Volume

Teichholz Formula

Teichholz LE, Kreulen T, Herman MV, Gorlin R. "Problems in Echocardiographic Volume Determinations: Echocardiographic-Angiographic Correlations in the Presence or Absence of Asynergy." *American Journal of Cardiology* 37(1):7-11, 1976.

Ejection Fraction

Pombo JF, Troy BL, Russell RO. "Left Ventricular Volumes and Ejection Fraction by Echocardiography." *Circulation* 43:480-490, 1971.

Fractional Shortening

Belenkie I, Nutter DO, Clark DW, McCraw DB, Raizner AE. "Assessment of Left Ventricular Dimensions and Function by Echocardiography." *American Journal of Cardiology* 31:755-762, June 1973.

Aorta/Left Atrium Ratio

Roelandt J. *Practical Echocardiography* in Ultrasound in Biomedicine, No. 1. Letchworth, Herts, England: Research Studies Press (a division of John Wiley & Sons Ltd), 1977, p. 270.

Schiller NB, Shah PM, Crawford M, et al. "Recommendations for Quantitation of the Left Ventricle by Two-Dimensional Echocardiography." *Journal of the American Society of Echocardiography* 2(5):364, 1989.

Heart Rate

Urdang, L. Ed. *Mosby's Medical & Nursing Dictionary*. St. Louis, Missouri: The C.V. Mosby Company, 1983, p. 492.

Stroke Index

Görge G, Erbel R, Brennecke R, Rupprecht HJ, Todt M, Meyer J. "High-Resolution Two-dimensional Echocardiography Improves the Quantification of Left Ventricular Function." *Journal of the American Society of Echocardiography* 5(2):125-134, 1992.

Stroke Volume

Görge G, Erbel R, Brennecke R, Rupprecht HJ, Todt M, Meyer J. "High-Resolution Two-dimensional Echocardiography Improves the Quantification of Left Ventricular Function." *Journal of the American Society of Echocardiography* 5(2):125-134, 1992.

Roelandt J. *Practical Echocardiography* in Ultrasound in Biomedicine, No. 1. Letchworth, Herts, England: Research Studies Press (a division of John Wiley & Sons Ltd), 1977, p. 270.

Doppler Mode Measurements and Calculations

Doppler References

Burns PN. "The Physical Principles of Doppler and Spectral Analysis." *Journal of Clinical Ultrasound* 15:567-590, 1987.

Acceleration Slope/Time

Gardin JM, Burn CS, Childs WJ, Henry WL. "Evaluation of blood flow velocity in the ascending aorta and main pulmonary artery of normal subjects by Doppler echocardiography." *American Heart Journal* 107(2):310-319, 1984.

Deceleration Slope/Time

Gardin JM, Burn CS, Childs WJ, Henry WL. "Evaluation of blood flow velocity in the ascending aorta and main pulmonary artery of normal subjects by Doppler echocardiography." *American Heart Journal* 107(2):310-319, 1984.

Pressure Half-Time

Hatle L, Angelsen B, Tromsdal A. "Noninvasive Assessment of Atrioventricular Pressure Half-time by Doppler Ultrasound." *Circulation* 60(5):1096-1103, 1979.

Mitral Valve Area Continuity Equation

Hatle L, Angelsen B, Tromsdal A. "Noninvasive Assessment of Atrioventricular Pressure Half-time by Doppler Ultrasound." *Circulation* 60(5):1096-1103, 1979.

Aortic Valve Area

Oh JK, Taliercio CP, Holmes DR, et. al. "Prediction of the Severity of Aortic Stenosis by Doppler Aortic Valve Area Determination: Prospective Doppler-Catheterization Correlation in 100 Patients." *Journal of the American College of Cardiology* 11(6):1227-1234, 1988.

Richards KL, Cannon SR, Miller JF, Crawford MH. "Calculation of aortic valve area by Doppler echocardiography: a direct application of the continuity equation." *Circulation* 73(5):964-969, 1986.

Heart Rate

Cohn PF, Brown EJ, Vlay SC. *Clinical Cardiovascular Physiology*. Philadelphia, Pennsylvania: W.B. Saunders Company, 1985, pp. 85-86.

Urdang, L. Ed. *Mosby's Medical & Nursing Dictionary*. St. Louis, Missouri: The C.V. Mosby Company, 1983, p. 492.

Cardiac Output

- Calafiore P, Stewart WJ. "Doppler Echocardiographic Quantitation of Volumetric Flow Rate." *Cardiology Clinics* 8(2):191-202, 1990.
- Haites NE, McLennan FM, Mowat DHR, Rawles JM. "Assessment of Cardiac Output by Ultrasound Technique Alone." *British Heart Journal* 53:123-129, 1985.
- Huntsman LL, Stewart DK, Barnes SR, Franklin SB, Colocousis JS, Hessel EA. "Noninvasive Doppler Determination of Cardiac Output in Man." *Circulation* 67(3):593-602, 1983.
- Ihlen H, Amlie JP, Dale J, et. al. "Determination of cardiac output by Doppler echocardiography." *British Heart Journal* 51:54-60, 1984.
- Magnin PA, Stewart JA, Myers S, von Ramm O, Kisslo JA. "Combined Doppler and Phased-array Echocardiographic Estimation of Cardiac Output." *Circulation* 63(2):388-392, 1981.

Mitral Valve, E point to A point Ratio

- Maron BJ, Spirito P, Green KJ, Wesley YE, Bonow RO, Arce J. "Noninvasive Assessment of Left Ventricular Diastolic Function by Pulsed Doppler Echocardiography in Patients with Hypertrophic Cardiomyopathy." *Journal of the American College of Cardiology* 10(4):733-742, 1987.

Mitral Valve, A point to E point Ratio

- Maron BJ, Spirito P, Green KJ, Wesley YE, Bonow RO, Arce J. "Noninvasive Assessment of Left Ventricular Diastolic Function by Pulsed Doppler Echocardiography in Patients with Hypertrophic Cardiomyopathy." *Journal of the American College of Cardiology* 10(4):733-742, 1987.

Mitral Valve, CA amplitude to CE amplitude Ratio

- Maron BJ, Spirito P, Green KJ, Wesley YE, Bonow RO, Arce J. "Noninvasive Assessment of Left Ventricular Diastolic Function by Pulsed Doppler Echocardiography in Patients with Hypertrophic Cardiomyopathy." *Journal of the American College of Cardiology* 10(4):733-742, 1987.

Isovolumic Relaxation Time

- Hatle L, Angelsen B, Tromsdal A. "Non-invasive estimation of pulmonary artery systolic pressure with Doppler ultrasound." *British Heart Journal* 45:157-165, 1981.
- Stevenson, JG. "Comparison of Several Noninvasive Methods for Estimation of Pulmonary Artery Pressure." *Journal of the American Society of Echocardiography* 2(3):157-166, 1989.

Left Ventricular Index of Myocardial Performance

- Tei C, Ling LH, Hodge DO, Bailey KR, Oh JK, Rodeheffer RJ, Tajik AJ, Seward JB. "New Index of Combined Systolic and Diastolic Myocardial Performance: A Simple and Reproducible Measure of Cardiac Function – A Study in Normals and Dilated Cardiomyopathy." *Journal of Cardiology* 26:357-366, 1995.

Right Ventricular Index of Myocardial Performance

- Tei C, Dujardin KS, Hodge DO, et al. "Doppler echocardiographic index for assessment of global right ventricular function." *Journal American Society of Echocardiography* 9:838-847, 1996.

Right Ventricular Systolic Pressure

- Yock PG, Popp RL. "Noninvasive estimation of right ventricular systolic pressure by Doppler ultrasound in patients with tricuspid regurgitation." *Circulation* 70(4):657-662, 1984.

Tricuspid Regurgitant Jet Plus Right Atrial Pressure

- Yock PG, Popp RL. "Noninvasive estimation of right ventricular systolic pressure by Doppler ultrasound in patients with tricuspid regurgitation." *Circulation* 70(4):657-662, 1984.

Maximum Pressure Gradient

Yoganathan AP, Cape EG, Sung H, Williams FP, Jimoh A. "Review of Hydrodynamic Principles for the Cardiologist: Applications to the Study of Blood Flow and Jets by Imaging Techniques." *Journal of the American College of Cardiology* 12(5):1344-1353, 1988.

Mean Pressure Gradient

Yoganathan AP, Cape EG, Sung H, Williams FP, Jimoh A. "Review of Hydrodynamic Principles for the Cardiologist: Applications to the Study of Blood Flow and Jets by Imaging Techniques." *Journal of the American College of Cardiology* 12(5):1344-1353, 1988.

Aortic Valve maximum Pressure Gradient

Yoganathan AP, Cape EG, Sung H, Williams FP, Jimoh A. "Review of Hydrodynamic Principles for the Cardiologist: Applications to the Study of Blood Flow and Jets by Imaging Techniques." *Journal of the American College of Cardiology* 12(5):1344-1353, 1988.

Aortic Valve mean Pressure Gradient

Yoganathan AP, Cape EG, Sung H, Williams FP, Jimoh A. "Review of Hydrodynamic Principles for the Cardiologist: Applications to the Study of Blood Flow and Jets by Imaging Techniques." *Journal of the American College of Cardiology* 12(5):1344-1353, 1988.

Mitral Valve maximum Pressure Gradient

Yoganathan AP, Cape EG, Sung H, Williams FP, Jimoh A. "Review of Hydrodynamic Principles for the Cardiologist: Applications to the Study of Blood Flow and Jets by Imaging Techniques." *Journal of the American College of Cardiology* 12(5):1344-1353, 1988.

Mitral Valve mean Pressure Gradient

Yoganathan AP, Cape EG, Sung H, Williams FP, Jimoh A. "Review of Hydrodynamic Principles for the Cardiologist: Applications to the Study of Blood Flow and Jets by Imaging Techniques." *Journal of the American College of Cardiology* 12(5):1344-1353, 1988.

Mitral Regurgitation maximum Pressure Gradient

Yoganathan AP, Cape EG, Sung H, Williams FP, Jimoh A. "Review of Hydrodynamic Principles for the Cardiologist: Applications to the Study of Blood Flow and Jets by Imaging Techniques." *Journal of the American College of Cardiology* 12(5):1344-1353, 1988.

Pulmonary Valve maximum Pressure Gradient

Yoganathan AP, Cape EG, Sung H, Williams FP, Jimoh A. "Review of Hydrodynamic Principles for the Cardiologist: Applications to the Study of Blood Flow and Jets by Imaging Techniques." *Journal of the American College of Cardiology* 12(5):1344-1353, 1988.

Pulmonary Valve mean Pressure Gradient

Yoganathan AP, Cape EG, Sung H, Williams FP, Jimoh A. "Review of Hydrodynamic Principles for the Cardiologist: Applications to the Study of Blood Flow and Jets by Imaging Techniques." *Journal of the American College of Cardiology* 12(5):1344-1353, 1988.

Tricuspid Regurgitation maximum Pressure Gradient

Yoganathan AP, Cape EG, Sung H, Williams FP, Jimoh A. "Review of Hydrodynamic Principles for the Cardiologist: Applications to the Study of Blood Flow and Jets by Imaging Techniques." *Journal of the American College of Cardiology* 12(5):1344-1353, 1988.

Tricuspid Regurgitation mean Pressure Gradient

Yoganathan AP, Cape EG, Sung H, Williams FP, Jimoh A. "Review of Hydrodynamic Principles for the Cardiologist: Applications to the Study of Blood Flow and Jets by Imaging Techniques." *Journal of the American College of Cardiology* 12(5):1344-1353, 1988.

Tricuspid Valve maximum Pressure Gradient

Yoganathan AP, Cape EG, Sung H, Williams FP, Jimoh A. "Review of Hydrodynamic Principles for the Cardiologist: Applications to the Study of Blood Flow and Jets by Imaging Techniques." *Journal of the American College of Cardiology* 12(5):1344-1353, 1988.

Tricuspid Valve mean Pressure Gradient

Yoganathan AP, Cape EG, Sung H, Williams FP, Jimoh A. "Review of Hydrodynamic Principles for the Cardiologist: Applications to the Study of Blood Flow and Jets by Imaging Techniques." *Journal of the American College of Cardiology* 12(5):1344-1353, 1988.

Left Ventricular Outflow Tract maximum Pressure Gradient

Yoganathan AP, Cape EG, Sung H, Williams FP, Jimoh A. "Review of Hydrodynamic Principles for the Cardiologist: Applications to the Study of Blood Flow and Jets by Imaging Techniques." *Journal of the American College of Cardiology* 12(5):1344-1353, 1988.

Left Ventricular Outflow Tract mean Pressure Gradient

Yoganathan AP, Cape EG, Sung H, Williams FP, Jimoh A. "Review of Hydrodynamic Principles for the Cardiologist: Applications to the Study of Blood Flow and Jets by Imaging Techniques." *Journal of the American College of Cardiology* 12(5):1344-1353, 1988.

Ventricular Septal Defect maximum Pressure Gradient

Yoganathan AP, Cape EG, Sung H, Williams FP, Jimoh A. "Review of Hydrodynamic Principles for the Cardiologist: Applications to the Study of Blood Flow and Jets by Imaging Techniques." *Journal of the American College of Cardiology* 12(5):1344-1353, 1988.

Additional Cardiac References

Left Atrial and Right Atrial Volumes

Bommer W, Weinert L, Neumann A, et al. "Determination of Right Atrial and Right Ventricular Size by Two-Dimensional Echocardiography." *Circulation* 60(1):91-100, 1979.

Lambertz H, Flachskampf FA, Heiliger R, et al. "New Echocardiographic and Angiographic Methods for Right Atrial Volume Determination: In Vitro Validation and In Vivo Results." *International Journal of Cardiac Imaging* 5(1):39-51, 1989.

Schabelman SE, Schiller NB, Anschultz RA, et al. "Comparison of Four Two-Dimensional Echocardiographic Views for Measuring Left Atrial Size." *American Journal of Cardiology* 41:391-6, 1978.

Wang Y, Gutman JM, Heilbron D, et al. "Atrial Volume in a Normal Adult Population by Two-Dimensional Echocardiography." *Chest* 86(4):595-601, 1984.

Left Ventricular Mass

- Bommer W, Weinert L, Neumann A, et al. "Determination of Right Atrial and Right Ventricular Size by Two-Dimensional Echocardiography." *Circulation* 60(1):91-100, 1979.
- Reichek N, Helak J, Plappert T, et al. "Anatomic Validation of Left Ventricular Mass Estimates from Clinical Two-Dimensional Echocardiography: Initial Results." *Circulation* 67(2):348-52, 1983.

Mitral Valve Area and Aortic Valve Area Using Continuity Equation with both the Peak Velocity and Velocity-time Integral

- Berger M, Berdoff RL, Gallerstein PE, et al. "Evaluation of Aortic Stenosis by Continuous Wave Doppler Ultrasound." *Journal of the American College of Cardiology* 3(1):150-6, 1984.
- Gossler KB, Goldberg SJ. "Velocity Gradients Across Normal Cardiac Valves." *American Journal of Cardiology* 67(1):99-101, 1991.
- Hatle L, Angelsen B. *Doppler Ultrasound in Cardiology: Physical Principles and Clinical Applications, 2nd edition*. Philadelphia, Pennsylvania: Lea & Febiger, 1985.
- Holen J, Aaslid R, Landmark K, et al. "Determination of Effective Orifice Area in Mitral Stenosis from Non-Invasive Ultrasound Doppler Data and Mitral Flow Rate." *Acta Medica Scandinavica* 201(1-2):83-8, 1977.
- Labovitz A. "Pressure Halftime, Deceleration Rate." *American Heart Journal* 8:1341-7, 1986.
- Loperfido F, Laurenzi F, Gimigliano F, et al. "A Comparison of the Assessment of Mitral Valve Area by Continuous Wave Doppler and by Cross Sectional Echocardiography." *British Heart Journal* 57(4):348-55, 1987.
- Nakatani S, Masuyama T, Kodama K, et al. "Value and Limitations of Doppler Echocardiography in the Quantification of Stenotic Mitral Valve Area: Comparison of the Pressure Half-Time and the Continuity Equation Methods." *Circulation* 77(1):78-85, 1988.
- Nishimura RA, Tajik AJ. "Quantitative Hemodynamics by Doppler Echocardiography: A Noninvasive Alternative to Cardiac Catheterization." *Progress in Cardiovascular Diseases* 36(4):309-42, 1994.
- Oh JK, Taliencio CP, Holmes DR, et al. "Prediction of the Severity of Aortic Stenosis by Doppler Aortic Valve Area Determination: Prospective Doppler-Catheterization Correlation in 100 Patients." *Journal of the American College of Cardiology* 11(6):1227-1234, 1988.
- Otto CM, Pearlman AS, Comess KA, et al. "Determination of the Stenotic Aortic Valve Area in Adults Using Doppler Echocardiography." *Journal of the American College of Cardiology* 7(3):509-17, 1986.
- Otto CM, Pearlman AS, Gardner CL, et al. "Simplification of the Doppler Continuity Equation for Calculating Stenotic Aortic Valve Area." *Journal of the American Society of Echocardiography* 1(2):155-7, 1988.
- Skjaerpe T, Hegrenaes L, Hatle L. "Noninvasive Estimation of Valve Area in Patients with Aortic Stenosis by Doppler Ultrasound and Two-Dimensional Echocardiography." *Circulation* 72(4):810-8, 1985.
- Stamm RB, Martin RP. "Quantification of Pressure Gradients Across Stenotic Valves by Doppler Ultrasound." *Journal of the American College of Cardiology* 2(4):707-18, 1983.
- Taylor R. "Evolution of the Continuity Equation in the Doppler Echocardiographic Assessment of the Severity of Valvular Aortic Stenosis." *Journal of the American Society of Echocardiography* 3(4):326-30, 1990.
- Thomas JD, Weyman AE. "Doppler Mitral Pressure Half-Time: A Clinical Tool in Search of Theoretical Justification." *Journal of the American College of Cardiology* 10(4):923-9, 1987.
- Wranne B, Ask P, Loyd D. "Analysis of Different Methods of Assessing the Stenotic Mitral Valve Area with Emphasis on the Pressure Gradient Half-Time Concept." *American Journal of Cardiology* 66(5):614-20, 1990.
- Young JB, Quinones MA, Waggoner AD, et al. "Diagnosis and Quantification of Aortic Stenosis with Pulsed Doppler Echocardiography." *American Journal of Cardiology* 45(5):987-94, 1980.
- Zoghbi WA, Farmer KL, Soto JG, et al. "Accurate Noninvasive Quantification of Stenotic Aortic Valve Area by Doppler Echocardiography." *Circulation* 73(3):452-9, 1986.

Cardiac Normal Range Values

- Baumgartner H, Hung J, Bermego J, Chambers JB, Evangelista A, Griffin BP, et al. "Echocardiographic Assessment of Valve Stenosis: EAE/ASE Recommendations for Clinical Practice." *Journal of the American Society of Echocardiography* 22(1):1-23, 2009.
- Lang RM, Bierig M, Devereux RB, Flachskampf FA, Foster E, Pellikka PA, et al. "Recommendations for Chamber Quantification." *Journal of the American Society of Echocardiography* 18(12):1440-1463, 2005.

Shortening / Wall Stress

- Colan SD, Borow KM, Newmann A. "Left ventricular end-systolic wall stress-velocity of fiber shortening relation: A load independent index of myocardial contractility." *Journal of the American College of Cardiology* 4:715-724, 1984.
- Engle SJ, DiSessa TG, Perloff JK, et al. "Mitral valve E point to ventricular septal separation in infants and children." *American Journal of Cardiology* 52:1084-1087, 1983.
- Ginzton LE, Kulick D. "Mitral valve E-point septal separation as an indicator of ejection fraction in patients with reversed septal motion." *Chest* 88:429-431, 1985.
- Massie BM, Schiller NB, Ratshin RA, et al. "Mitral-septal separation: A new echocardiographic index of left ventricular function." *American Journal of Cardiology* 39:1008-1016, 1977.
- Quinones MA, Mokotoff DM, Nouri S, et al. "Noninvasive quantification of left ventricular wall stress: Validation of method and application to assessment of chronic pressure overload." *American Journal of Cardiology* 45:782-790, 1980.
- Reichek N, Wilson J, St. John Sutton M, Plappert TA, Goldberg S, Hirshfeld J. "Noninvasive Determination of Left Ventricular End-systolic Stress: Validation of the Method and Initial Application." *Circulation* 65(1), 1982.
- Sahn DJ, Deely WJ, Hagen AD, et al. "Echocardiographic assessment of left ventricular performance in normal newborns." *Circulation* 49:232-236, 1974.

Relative Wall Thickness

- Lang RM, et al. "Recommendations for chamber quantification." *Eur J Echocardiography* 7(2): 79-108, 2006.

syngo VVI Clinical References

- Butz T, Lang CN, van Bracht M, Prull MW, Yeni H, Maagh P, Plehn G, Meissner A, Trappe HJ. "Segment-orientated Analysis of Two-dimensional Strain and Strain Rate as Assessed by Velocity Vector Imaging in Patients with Acute Myocardial Infarction." *International Journal of Medical Sciences* 8(2):106–113, 2011.
- Pirat B, Khoury DS, Hartley CJ, Tiller L, Rao L, Schulz DG, Nagueh SF, Zoghbi WA. "A Novel Feature-tracking Echocardiographic Method for the Quantitation of Regional Myocardial Function." *Journal of the American College of Cardiology* 51(6):651–659, 2008.
- Rodríguez-Bailón I, Jiménez-Navarro MF, Pérez-González R, García-Orta R, Morillo-Velarde E, de Teresa-Galvána E. "Left Ventricular Deformation and Two-Dimensional Echocardiography: Temporal and Other Parameter Values in Normal Subjects." *Revista Española de Cardiología* 63(10):1195–1199, 2010.
- Willruth AM, Geipel A, Berg C, Fimmers R, Gembruch U. "Assessment of Left Ventricular Global and Regional Longitudinal Peak Systolic Strain, Strain Rate and Velocity with Feature Tracking in Healthy Fetuses." *Ultraschall in der Medizin* 33(7):E293–298, 2012. Epub 2011 Feb 3.

Appendix H Data File Descriptions

- syngo Velocity Vector Imaging..... 2**
 - Bullseye Summary Results..... 2
 - Time to Peak Summary Results 2
 - Velocity Results 2
 - Displacement Results 3
 - Strain Results 4
 - Strain Rate Results 5
 - Volume and Area Results 6
 - Descriptions for Exported Data Files 7

syngo Velocity Vector Imaging

See also: For information about syngo Velocity Vector Imaging, refer to Chapter D4 in the Clinical Applications Guide.

Bullseye Summary Results

The displayed results depend on the selected view, parameter, layer, and direction of analysis.

Time to Peak Summary Results

Label	Description
Maximum Opposing Wall Delay	Maximum opposing wall delay (MOWD)
Maximum Wall Delay	(Available only in SAX views) Maximum opposing wall delay (MOWD)
Average	Average of all segments

Velocity Results

Result Label	Description	Units
Pk cm/s	Peak systolic velocity	cm/s
Pk deg/s	(Available only in SAX views) Peak systolic velocity	deg/s
ES	End systolic velocity	%
TTP ms	Time-to-peak systolic velocity	ms
es-gcs endo	(Available only in Apical views or SAX views) End systolic global circumferential strain endocardium	%
es-gcs-myocardium	(Available only in Apical views or SAX views when the Endo+Epi selection is enabled) End systolic global circumferential strain myocardium	%
es-rot endo	(Available only in SAX views) End systolic rotational endocardium	deg
es-rot-myocardium	(Available only in SAX views when the Endo+Epi selection is enabled) End systolic rotational myocardium	deg
es-gls	(Available only in atrium) End systolic global longitudinal strain	%
es-gls endo	(Available only in Apical views or right ventricle) End systolic global longitudinal endocardium	%
es-gls myocardium	(Available only in Apical views or right ventricle when the Endo+Epi selection is enabled) End systolic global longitudinal myocardium	%

Displacement Results

Result Label	Description	Units
Pk mm	Peak displacement	mm
Pk deg	(Available only in SAX views) Peak displacement	deg
ES	End systolic displacement	%
TTP ms	Time-to-peak displacement	ms
es-gcs endo	(Available only in Apical views or SAX views) End systolic global circumferential strain endocardium	%
es-gcs-myo	(Available only in Apical views or SAX views when the Endo+Epi selection is enabled) End systolic global circumferential strain myocardium	%
es-rot endo	(Available only in SAX views) End systolic rotational endocardium	deg
es-rot-myo	(Available only in SAX views when the Endo+Epi selection is enabled) End systolic rotational myocardium	deg
es-gls	(Available only in atrium) End systolic global longitudinal strain	%
es-gls endo	(Available only in Apical views or right ventricle) End systolic global longitudinal endocardium	%
es-gls myo	(Available only in Apical views or right ventricle when the Endo+Epi selection is enabled) End systolic global longitudinal myocardium	%

Strain Results

Result Label	Description	Units
Pk %	Peak systolic strain	%
ES	End systolic strain	%
TTP ms	Time-to-peak strain	ms
es-gls	(Available only in atrium) End systolic global longitudinal strain	%
es-gls endo	(Available only in Apical views or right ventricle) End systolic global longitudinal strain endocardium	%
es-gls myo	(Available only in Apical views or right ventricle when the Endo+Epi selection is enabled) End systolic global longitudinal strain myocardium	%
es-gcs endo	(Available only in Apical views or SAX views) End systolic global circumferential strain endocardium	%
es-gcs-myocardium	(Available only in Apical views or SAX views when the Endo+Epi selection is enabled) End systolic global circumferential strain myocardium	%
es-rot endo	(Available only in SAX views) End systolic rotational endocardium	deg
es-rot-myocardium	(Available only in SAX views when the Endo+Epi selection is enabled) End systolic rotational myocardium	deg

Strain Rate Results

Result Label	Description	Units
Pk 1/s	Peak strain rate	1/s
ES	End systolic strain rate	%
TTP ms	Time-to-peak systolic strain rate	ms
es-gcs endo	(Available only in Apical views or SAX views) End systolic global circumferential strain endocardium	%
es-gcs-myo	(Available only Apical views or SAX views when the Endo+Epi selection is enabled) End systolic global circumferential strain myocardium	%
es-rot endo	(Available only in SAX views) End systolic rotational endocardium	deg
es-rot-myo	(Available only in SAX views when the Endo+Epi selection is enabled) End systolic rotational myocardium	deg
es-gls	(Available only in atrium) End systolic global longitudinal strain	%
es-gls endo	(Available only in Apical views or right ventricle) End systolic global longitudinal endocardium	%
es-gls myo	(Available only in Apical views or right ventricle when the Endo+Epi selection is enabled) End systolic global longitudinal myocardium	%

Volume and Area Results

Result Label	Description	Units
EF	Ejection fraction	%
EDV	End-diastolic volume	ml
ESV	End-systolic volume	ml
FAC	Fractional area change	%
MyoGCS	Myocardium global circumferential strain	%
MyoGLS	Myocardium global longitudinal strain	%
EndoGCS	Endocardium global circumferential strain	%
EndoGLS	Endocardium global longitudinal strain	%
EDA	End diastolic area	cm ²
ESA	End systolic area	cm ²
EndoRot	Endocardial rotation	deg
MyoRot	Myocardium rotation	deg
GRS	Global radial strain	%
SD-RS-Peak	Standard deviation of peak radial strain	%
SD-RS-Syst.	Standard deviation of radial strain at systole	%
SD-CS-Peak	Standard deviation of peak circumferential strain	%
SD-CS-Syst.	Standard deviation of circumferential strain at systole	%
SD-Ph%-RS	Standard deviation of phase of radial strain	%
SD-Ph%-CS	Standard deviation of phase of circumferential strain	%
SD-Ph%-TS	Standard deviation of phase of transverse strain	%
SD-Ph%-LS	Standard deviation of phase of longitudinal strain	%
SD-TS-Peak	Standard deviation of peak transverse strain	%
SD-TS-Syst.	Standard deviation of transverse strain at systole	%
SD-LS-Peak	Standard deviation of peak longitudinal strain	%
SD-LS-Syst.	Standard deviation of longitudinal strain at systole	%
SD-Ttp%-TS	Standard deviation of time to peak of transverse strain	%
SD-Ttp%-LS	Standard deviation of time to peak of longitudinal strain	%
SD-Ttp%-RS	Standard deviation of time to peak of radial strain	%
SD-Ttp%-CS	Standard deviation of time to peak of circumferential strain	%
Delta-ROT	Delta of rotational strain	deg
FreeWallSt	Free wall strain	%
EDL	End-diastolic length	cm
ESL	End-systolic length	cm
EDDbas	End-diastolic diameter basal	cm
ESDbas	End-systolic diameter basal	cm

Descriptions for Exported Data Files

The program stores exported files in an external USB-compatible storage device.

Selection	Processing Algorithm(s)	Description
Axis	All	Axis for the selected processing algorithm, for example, Long or Short.
View	All	Selected view, for example, right ventricle.
segments model	All	Segment information, for example, apical segment.
eD	All	End diastole.
eS	All	End systole.
Bpm	All	Heart rate, in beats per minute.
PixelDimension	All	Millimeters per pixel.
TracePointsEpiXY (pixel)	Apical for Endo/Epi, Short Axis for Endo/Epi	(x1,y1,x2,y2,x3,y3,...) pixel coordinates of the endo and epi user-defined points to be tracked.
TracePointsEpiXY (mm)	Apical for Endo/Epi, Short Axis for Endo/Epi	(x1,y1,x2,y2,x3,y3,...) coordinates of the endo and epi user-defined points to be tracked, converted to millimeters.
TracePointsEndoXY (pixel)	Apical for Endo/Epi, Short Axis for Endo/Epi	(x1,y1,x2,y2,x3,y3,...) pixel coordinates of the endo and epi user-defined points to be tracked.
TracePointsEndoXY (mm)	Apical for Endo/Epi, Short Axis for Endo/Epi	(x1,y1,x2,y2,x3,y3,...) coordinates of the endo and epi user-defined points to be tracked, converted to millimeters.
EcgPoints	All	Amplitude of the electrocardiogram, sampled every two milliseconds.
RWaves	All	Amplitude of the R peak from QRS complex.
FrameTime (msec)	All	Time between adjacent frames.
TimeProgression (msec)	All	Time since the first frame.
Time (msec)	All	Time for the selected location.
Velocity (cm/sec)	All	Velocity towards (positive) or away from (negative) the velocity reference point.

Selection	Processing Algorithm(s)	Description
Strain (%)	All	Strain (percentage) for all points along the trace, for every frame.
StrainRate (1/s)	All	Strain rate (per second) for all points along the trace, for every frame.
EndoX (pixel coord)	All	x-coordinates of all endo points along the trace, for every frame.
EndoY (pixel coord)	All	y-coordinates of all endo points along the trace, for every frame.
EndoVelX (cm/sec)	All	x component of the velocity for every endo point along the trace, for every frame.
EndoVelY (cm/sec)	All	y component of the velocity for every endo point along the trace, for every frame.
EpiX (pixel coord)	Apical for Endo/Epi, Short Axis for Endo/Epi	x-coordinates of the epi points along the trace, for every frame.
EpiY (pixel coord)	Apical for Endo/Epi, Short Axis for Endo/Epi	y-coordinates of the epi points along the trace, for every frame.
EpiVelX (cm/sec)	Apical for Endo/Epi, Short Axis for Endo/Epi	x component of the velocity for epi points along the trace, for every frame.
EpiVelY (cm/sec)	Apical for Endo/Epi, Short Axis for Endo/Epi	y component of the velocity for epi points along the trace, for every frame.
RadialVelocity (cm/s)	Short Axis, Short Axis for Endo/Epi	Radial velocity along the trace, for every frame.
RotationVelocity (deg/s)	Short Axis, Short Axis for Endo/Epi	Rotation velocity along the trace, for every frame.
RadialDisplacement (mm)	Short Axis, Short Axis for Endo/Epi	Radial displacement along the trace, for every frame.
RotationDisplacement (deg)	Short Axis, Short Axis for Endo/Epi	Rotation displacement along the trace, for every frame.
RadialStrain (%)	Short Axis for Endo/Epi	Radial strain (percentage) along the trace, for every frame.
CircumferentialStrain (%)	Short Axis, Short Axis for Endo/Epi	Circumferential strain (percentage) along the trace, for every frame.
RadialStrainRate (1/s)	Short Axis for Endo/Epi	Radial strain rate (per second) along the trace, for every frame.
CircumferentialStrainRate (1/s)	Short Axis, Short Axis for Endo/Epi	Circumferential strain rate (per second) along the trace, for every frame.

Selection	Processing Algorithm(s)	Description
Shear (%)	Apical for Endo/Epi, Short Axis for Endo/Epi, Right Ventricle for Endo/Epi	The difference in the tangential displacement of the endo and epi points (percentage), divided by the endo-epi distance.
ShearRate (1/s)	Apical for Endo/Epi, Short Axis for Endo/Epi, Right Ventricle for Endo/Epi	The difference in the tangential velocity of the endo and epi points (per second), divided by the endo-epi distance.
Vol (ml)	Apical, Apical for Endo/Epi	Modified method of disk volume for every frame, converted to milliliters.
Area (cm ²)	Short Axis, Short Axis for Endo/Epi, Right Ventricle for Endo/Epi	Area for every frame, in squared centimeters.
dV/dt (ml/s)	Apical, Apical for Endo/Epi	The rate of change in the volume.
EndoEpiDistance (mm)	Apical for Endo/Epi, Short Axis for Endo/Epi, Right Ventricle for Endo/Epi	Distance between corresponding endo and epi points. Not the perpendicular distance between the point and the other trace.
EndoGLS (%)	Apical, Apical for Endo/Epi, Right Ventricle for Endo/Epi, Atrium	Endocardium global longitudinal strain (percentage) for the endo points along the trace, for every frame.
EndoGCS (%)	Apical, Apical for Endo/Epi, Short Axis, Short Axis for Endo/Epi	Endocardium global circumferential strain (percentage) for the endo points along the trace, for every frame.
MyoGLS (%)	Apical for Endo/Epi, Right Ventricle for Endo/Epi	Myocardium global longitudinal strain (percentage) for the endo and epi points along the trace, for every frame.
MyoGCS (%)	Apical for Endo/Epi, Short Axis for Endo/Epi	Myocardium global circumferential strain (percentage) for the endo and epi points along the trace, for every frame.
GRS (%)	Apical for Endo/Epi, Short Axis for Endo/Epi, Right Ventricle for Endo/Epi	Global radial strain (percentage) for the endo and epi points along the trace, for every frame.
TransverseVelocity (cm/s)	Apical, Apical for Endo/Epi, Atrium, Right Ventricle for Endo/Epi	Transverse velocity for the endo and epi points along the trace, for every frame.

Selection	Processing Algorithm(s)	Description
LongitudinalVelocity (cm/s)	Apical, Apical for Endo/Epi, Atrium, Right Ventricle for Endo/Epi	Longitudinal velocity for the endo and epi points along the trace, for every frame.
TransverseDisplacement (mm)	Apical, Apical for Endo/Epi, Atrium, Right Ventricle for Endo/Epi	Transverse displacement for the endo and epi points along the trace, for every frame.
LongitudinalDisplacement (mm)	Apical, Apical for Endo/Epi, Atrium, Right Ventricle for Endo/Epi	Longitudinal displacement for the endo and epi points along the trace, for every frame.
TransverseStrain (%)	Apical for Endo/Epi, Right Ventricle for Endo/Epi	Transverse strain (percentage) for the endo and epi points along the trace, for every frame.
LongitudinalStrain (%)	Apical, Apical for Endo/Epi, Atrium, Right Ventricle for Endo/Epi	Longitudinal strain (percentage) for the endo and epi points along the trace, for every frame.
TransverseStrainRate (1/s)	Apical for Endo/Epi, Right Ventricle for Endo/Epi	Transverse strain rate (per second) for the endo and epi points along the trace, for every frame.
LongitudinalStrainRate (1/s)	Apical, Apical for Endo/Epi, Atrium, Right Ventricle for Endo/Epi	Longitudinal strain rate (per second) for the endo and epi points along the trace, for every frame.
EndoRot (deg)	Short Axis, Short Axis for Endo/Epi	Endo rotational degree.
MyoRot (deg)	Short Axis for Endo/Epi	Myocardium rotational degree.

