XEMA

XEMA CO., Ltd.
bldg. 48/4, 9th Parkovaya str., 105264, Moscow, Russia
Tel./Fax: +7 (495) 510-57-07, +7 (495) 737-39-36
E-mail: info@xema.ru, info@xema-medica.com
Internet: www.xema.ru, www.xema-medica.com

STATEMENT

We, XEMA Co., Ltd. having a registered office at 48, 9th Parkovaya st, 104264 Moscow, Russia, assign Sanmedico Srl. having a registered office at str. A. Corobceanu 7A, apt. 9, Chişinău MD 2012, Moldova, as authorized representative in correspondence with the conditions of directive 93/42/EEC, 98/79/EEC and 90/385/EEC .

We declare that the company mentioned above is authorized to register, notify, renew or modify the registration of medical devices on the territory of the Republic of Moldova.

Date : November, 29, 2017

Signature:

Andrei P. Redkin

Deputy general manager

Certificate

Of Marketing Authorization of Medical Product

Nr. **AR/IVMD/Xema/01/2020**

Issued on the basis of the Declaration of conformity and registration taking into account Article 10 of Directive 98/79/EC on In Vitro Diagnostic Medical Devices and Medical Devices Act (MPG) § § 5,25,29,30

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Valid to : 2022-05-25
Gültig bis

Represented in the EC by Polmed.de
Steinacker 5, 73773 Aichwald, Germany
email: info@polmed.de
tel: +49 711 52853279

Polmed.de

Annex to the Certificate No.:

Anhang zum Zertifikat Nr.:

AR/IVMD/Xema/01/2020

The following medical devices can be placed on the market in the Federal Republic of Germany, in the member states of the European Economic Community (EEC) and in the other contract states of the agreement about the European Economic Area.

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
Represented in the EC by Polmed.de
 Steinacker 5, 73773 Aichwald, Germany
 email: info@polmed.de
 tel: +49 711 52853279

Date: January 01, 2020

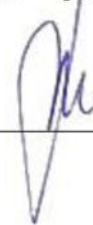Polmed.de

MANAGEMENT SYSTEM CERTIFICATE

Certificate No:
282710-2019-AQ-MCW-FINAS

Initial certification date:
14 February 2019

Valid:
14 February 2019 - 14 February 2022

This is to certify that the management system of

XEMA Co, LTD

bldg. 48, 9-th Parkovaya str., Moscow, Russian Federation, 105264
and the sites as mentioned in the appendix accompanying this certificate

has been found to conform to the Quality Management System standard:
ISO 9001:2015

This certificate is valid for the following scope:

Design and development, manufacturing and sales of in vitro tests for food and feed control, clinical and veterinary diagnostics and forensic investigations.

Place and date:
Moscow, 14 February 2019

For the issuing office:
DNV GL – Business Assurance
Trekhpudny per. 9 build. 2, office 406,
Moscow, Russian Federation

S. Groubine

Serguei Groubine
Management Representative

Certificate No: 282710-2019-AQ-MCW-FINAS
Place and date: Moscow, 14 February 2019

Appendix to Certificate

XEMA Co, LTD

Locations included in the certification are as follows:

Site Name	Site Address	Site Scope
XEMA Co, LTD	bldg. 48, 9-th Parkovaya str., Moscow, Russian Federation, 105264	Design and development, manufacturing and sales of in vitro tests for food and feed control, clinical and veterinary diagnostics and forensic investigations.
XEMA Co, LTD (Production site)	Trubetskaya str., 2B, Balashikha, Moscow region, Russian Federation, 125000	Design and development, manufacturing and sales of in vitro tests for food and feed control, clinical and veterinary diagnostics and forensic investigations.

MANAGEMENT SYSTEM CERTIFICATE

Certificate No:
53899-2009-AQ-MCW-FINAS

Initial certification date:
22 May 2009

Valid:
22 January 2019 - 30 April 2021

This is to certify that the management system of

XEMA CO., LTD.

bldg. 48, 9-th Parkovaya str., Moscow, Russian Federation, 105264
and the sites as mentioned in the appendix accompanying this certificate

has been found to conform to the Quality Management System standard:
ISO 13485:2016

This certificate is valid for the following scope:

**DESIGN, DEVELOPMENT, MANUFACTURING AND SALES OF KITS FOR IVD
USE.**

Place and date:
Moscow, 22 January 2019

For the issuing office:
DNV GL – Business Assurance
Trekhpudny per. 9 build. 2, office 406,
Moscow, Russian Federation

S. Groubine

Serguei Groubine
Management Representative

Certificate No: 53899-2009-AQ-MCW-FINAS
Place and date: Moscow, 22 January 2019

Appendix to Certificate

XEMA CO., LTD.

Locations included in the certification are as follows:

Site Name	Site Address	Site Scope
XEMA CO., LTD.	bldg. 48, 9-th Parkovaya str., Moscow, Russian Federation, 105264	DESIGN, DEVELOPMENT, MANUFACTURING AND SALES OF KITS FOR IVD USE.
XEMA Co., LTD (production site)	Trubetskaya str., 2B, Balashikha, Moscow region, Russian Federation, 125000	DESIGN, DEVELOPMENT, MANUFACTURING AND SALES OF KITS FOR IVD USE.

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**A SOLID-PHASE ENZYME IMMUNOASSAY
FOR THE QUANTITATIVE DETERMINATION
OF AUTOANTIBODIES TO THYROGLOBULIN
IN HUMAN SERUM OR PLASMA**

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
Authorized Representative in EU:
Polmed.de
Steinacker 20, D-73773
Aichwald, Germany
e-mail: info@polmed.de

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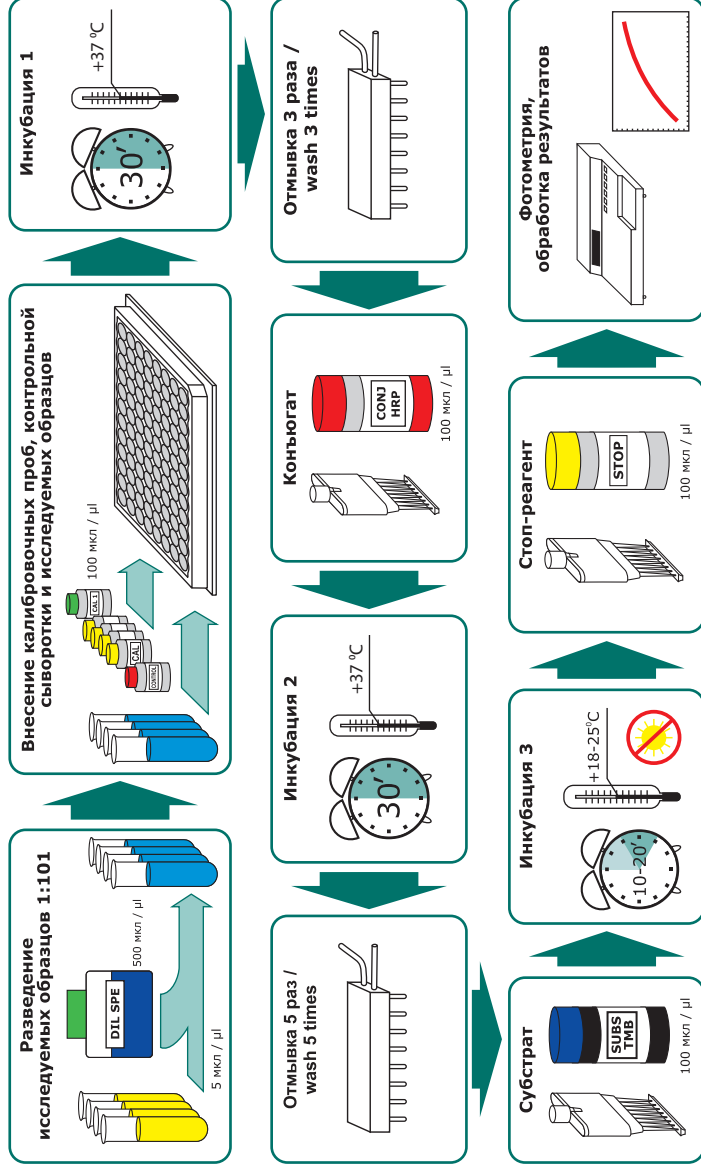CONTENT

1. INTENDED USE	10
2. SUMMARY AND EXPLANATION	10
3. PRINCIPLE OF THE TEST	10
4. WARNINGS AND PRECAUTIONS	11
5. KIT COMPONENTS	12
6. SPECIMEN COLLECTION AND STORAGE	13
7. TEST PROCEDURE	13
8. QUALITY CONTROL	15
9. CALCULATION OF RESULTS	15
10. EXPECTED VALUES	15
11. PERFORMANCE CHARACTERISTICS	16
12. LITERATURE	16

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Instruction for use

A SOLID-PHASE ENZYME IMMUNOASSAY FOR THE QUANTITATIVE DETERMINATION OF AUTOANTIBODIES TO THYROGLOBULIN IN HUMAN SERUM OR PLASMA

1. INTENDED USE

A solid-phase enzyme immunoassay for the quantitative determination of autoantibodies to thyroglobulin in serum or plasma.

This kit is designed for measurement of autoantibodies to thyroglobulin in serum or plasma. For possibility of use with other sample types, please, refer to Application Notes (on request). The kit contains reagents sufficient for 96 determinations and allows to analyze 42 unknown samples in duplicates.

2. SUMMARY AND EXPLANATION

Thyroglobulin (TG) is a well known target for autoantibodies occurring in thyroid autoimmunity (Graves' disease and Hashimoto's thyroiditis). Anti-TG antibodies mostly belong to the IgG class. Low to moderate levels of anti-TG antibodies can be found in sera of other autoimmune patients (eg systemic lupus erythematosus or Sjogren syndrome). In some cases anti-TG positive sera may show negativity for other type of anti-thyroid antibodies – anti-TPO. Therefore, combined determination of both types of anti-thyroid antibodies (anti-TPO + anti-TG) provides most sensitive laboratory diagnostic tool for thyroid autoimmunity. Separately from autoimmunity, anti-TG antibodies may develop in patients suffering from thyroid cancer. High level of anti-TG in such patients may interfere with correct determination of serum thyroglobulin which serves as tumour marker for therapy control in this group of patients.

3. PRINCIPLE OF THE TEST

This test is based on indirect enzyme immunoassay principle. Tested specimen is placed into the microwells coated by the antigen. Antibodies from the specimen bind coated antigen on the microwell surface. Unbound material is removed by washing procedure. Second antibodies directed towards species specific Ig, labelled with peroxidase enzyme, are then added into the microwells. After subsequent washing procedure, the remaining enzymatic activity bound to the microwell surface is detected and quantified by addition of chromogen-substrate mixture, stop solution and photometry at 450 nm. Optical density in the microwell is directly related to the quantity of the measured analyte in the specimen.

4. WARNINGS AND PRECAUTIONS

4.1. For professional use only.

4.2. This kit is intended for in vitro diagnostic use only.

4.3. INFECTION HAZARD: There is no available test methods that can absolutely assure that Hepatitis B and C viruses, HIV-1/2, or other infectious agents are not present in the reagents of this kit. All human products, including patient samples, should be considered potentially infectious. Handling and disposal should be in accordance with the procedures defined by an appropriate national biohazard safety guidelines or regulations.

4.4. Avoid contact with stop solution containing 5.0% H_2SO_4 . It may cause skin irritation and burns.

4.5. Wear disposable latex gloves when handling specimens and reagents. Microbial contamination of reagents may give false results.

4.6. Do not use the kit beyond the expiration date.

4.7. All indicated volumes have to be performed according to the protocol. Optimal test results are only obtained when using calibrated pipettes and microplate readers.

4.8. Do not smoke, eat, drink or apply cosmetics in areas where specimens or kit reagents are handled.

4.9. Chemicals and prepared or used reagents have to be treated as hazardous waste according to the national biohazard safety guidelines or regulations.

4.10. Do not mix reagents from different lots.

4.11. Replace caps on reagents immediately. Do not swap caps.

4.12. Do not pipette reagents by mouth.

4.13. Specimens must not contain any AZIDE compounds – they inhibit activity of peroxidase.

4.14. Material Safety Data Sheet for this product is available upon request directly from XEMA Co., Ltd.

4.15. The Material Safety Data Sheet fit the requirements of EU Guideline 91/155 EC.

5.1. Contents of the Kit

5. KIT COMPONENTS

Symbol	Description	Qty	Units	Colour code	Stability of opened/diluted components
1 SORB MTP	aTG EIA strips, 8x12 wells	1	pcs		until exp.date
2 CAL 1-5	Calibrator set, 1.1 ml each. The set contains 5 calibrators: 0; 100; 300; 1000; 3000 IU/ml	5	pcs	blue (C1 – colourless)	2 months
3 CONTROL	Control serum (1.1 ml)	1	pcs	colourless	2 months
4 CONJ HRP	Conjugate, 14 ml	1	pcs	purple	until exp.date
5 DIL SPE	EIA buffer 50 ml	1	pcs	blue	until exp.date
6 SUBS TMB	Substrate solution, 14 ml	1	pcs	colourless	until exp.date
7 BUF WASH 26X	Washing solution concentrate 26X, 22 ml	1	pcs	colourless	Concentrate – until exp.date Diluted washing solution – 45 days at 2-8 °C or 15 days at RT
8 STOP	Stop solution, 14 ml	1	pcs	colourless	until exp.date
9 N003	Plate sealing tape	2	pcs		N/A
10 K132I	Instruction aTG EIA	1	pcs		N/A
11 K132Q	QC data sheet aTG EIA	1	pcs		N/A

5.2. Equipment and material required but not provided

- Distilled or deionized water;
- Automatic or semiautomatic multichannel micropipettes, 100–250 µl, is useful but not essential;
- Calibrated micropipettes with variable volume, range volume 5–250 µl;
- Calibrated microplate photometer with 450 nm wavelength and OD measuring range 0–3.0.

5.3. Storage and stability of the Kit

Store the whole kit at +2...+8 °C upon receipt until the expiration date.

After opening the pouch keep unused microtiter wells **TIGHTLY SEALED BY ADHESIVE TAPE (INCLUDED)** to minimize exposure to moisture.

6. SPECIMEN COLLECTION AND STORAGE

This kit is intended for use with serum or plasma (ACD- or heparinized). Grossly hemolytic, lipemic, or turbid samples should be avoided.

Specimens may be stored for up to 48 hours at +2...+8 °C before testing. For a longer storage, the specimens should be frozen at -20 °C or lower. Repeated freezing/thawing should be avoided.

7. TEST PROCEDURE**7.1. Reagent Preparation**

- All reagents (including unsealed microstrips) should be allowed to reach room temperature (+18...+25 °C) before use.
- All reagents should be mixed by gentle inversion or vortexing prior to use. Avoid foam formation.
- It is recommended to spin down shortly the tubes with calibrators on low speed centrifuge.
- Prepare washing solution from the concentrate BUF WASH 26X by 26 dilutions in distilled water.

7.2. Procedural Note:

It is recommended that pipetting of all calibrators and samples should be completed within 3 minutes.

7.3. Assay flowchart

See the example of calibration graphic in Quality Control data sheet.

7.4. Assay procedure

1	Put the desired number of microstrips into the frame; allocate 12 wells for the calibrators CAL 1-5 and control samples CONTROL and two wells for each unknown sample. DO NOT REMOVE ADHESIVE SEALING TAPE FROM UNUSED STRIPS.
2	Dilute all samples using buffer DIL SPE (EIA buffer) 101 fold. (5 µl of sample + 500 µl of diluent) . Do not dilute control sample and calibrators.
3	If suggested analyte concentration in the sample exceeds the highest calibrator, additionally dilute this sample accordingly, using DIL SPE (EIA buffer). Use of other buffers or reagents for sample dilution may lead to incorrect measurement.
4	Pipet 100 µl of calibrators CAL 1-5, control sample CONTROL and unknown samples into the wells. Cover the wells by plate adhesive tape (included into the kit).
5	Incubate 30 minutes at +37 °C.
6	Prepare washing solution by 26X dilution of washing solution concentrate BUF WASH 26X by distilled water. Minimal quantity of washing solution should be 250 µl per well. Wash strips 3 times.
7	Dispense 100 µl of CONJ HRP into the wells. Cover the wells by plate adhesive tape.
8	Incubate 30 minutes at +37 °C.
9	Wash the strips 5 times.
10	Dispense 100 µl of SUBS TMB into the wells.
11	Incubate 10-20 minutes at +18...+25 °C.
12	Dispense 100 µl of STOP into the wells.
13	Measure OD (optical density) at 450 nm.
14	Set photometer blank on first calibrator.
15	Apply point-by-point method for data reduction.

7.5. Sample processing

Material type	Notes on material collection, storage and handling	Sample dilution example	Sample into the well, µl	Calculation factor
serum or plasma	Grossly hemolytic, lipemic, or turbid samples should be avoided and should be treated by centrifugation before testing.	5 µl of sample + 500 µl of diluent	100	1

8. QUALITY CONTROL

It is recommended to use control samples according to state and federal regulations. The use of control samples is advised to assure the day to day validity of results.

The test must be performed exactly as per the manufacturer's instructions for use. Moreover the user must strictly adhere to the rules of GLP (Good Laboratory Practice) or other applicable federal, state, and local standards and/or laws. This is especially relevant for the use of control reagents. It is important to always include, within the test procedure, a sufficient number of controls for validating the accuracy and precision of the test.

The test results are valid only if all controls are within the specified ranges and if all other test parameters are also within the given assay specifications.

9. CALCULATION OF RESULTS


9.1. Calculate the mean absorbance values (OD450) for each pair of calibrators and samples.

9.2. Plot a calibration curve on graph paper: OD versus autoantibodies to thyroglobulin concentration.

9.3. Determine the corresponding concentration of autoantibodies to thyroglobulin in unknown samples from the calibration curve. Manual or computerized data reduction is applicable on this stage. Point-by-point or linear data reduction is recommended due to non-linear shape of curve.

9.4. Below is presented a typical example of a standard curve with the XEMA Co. Not for calculations!

Calibrators	Value	Absorbance Units (450 nm)
CAL 1	0 IU/ml	0.04
CAL 2	100 IU/ml	0.15
CAL 3	300 IU/ml	0.35
CAL 4	1000 IU/ml	1.07
CAL 5	3000 IU/ml	2.30

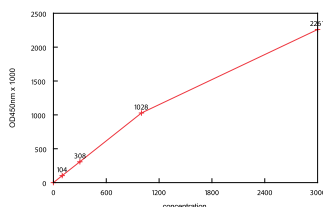10. EXPECTED VALUES

Therapeutical consequences should not be based on results of IVD methods alone – all available clinical and laboratory findings should be used by a physician to elaborate therapeutically measures. Each laboratory should establish its own normal range for aTG. Based on data obtained by XEMA, the following normal range is recommended (see below).

Sex, age	Units, IU/ml	
	Lower limit	Upper limit
Males	-	100
Females	-	100
>50 yrs	-	150

11. PERFORMANCE CHARACTERISTICS

11.1. Analytical sensitivity

Sensitivity of the assay was assessed as being 5.0 IU/ml.

11.2. Linearity

Linearity was checked by assaying dilution series of 5 samples with different autoantibodies to thyroglobulin concentrations. Linearity percentages obtained ranged within 90 to 110%.

11.3. Recovery

Recovery was estimated by assaying 5 mixed samples with known autoantibodies to thyroglobulin concentrations. The recovery percentages ranged from 90 to 110%.

12. LITERATURE

1. U Feldt-Rasmussen – Analytical and clinical performance goals for testing autoantibodies to thyroperoxidase, thyroglobulin, and thyrotropin receptor. Clin. Chem., Jan 1996; 42: 160 – 163.
2. PW Ladenson – Optimal laboratory testing for diagnosis and monitoring of thyroid nodules, goiter, and thyroid cancer. Clin. Chem., Jan 1996; 42: 183 – 187.
3. Anthony P. Weetman – Graves' Disease. N. Engl. J. Med., Oct 2000; 343: 1236 – 1248.

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