



Viale G. Di Vittorio, 30
20048 Pantigliate (MI) Italy
Phone. +39 02.55230061
www.histoline.com
email: export@histoline.com

GEMINI Semi-Automatic Rotary Microtome Data Sheet



Advancing Science & Medicine Through Histo-Cytology



Manufacturer: Histo-Line Laboratories s.r.l

In vitro diagnostic -medical device

Version 1.0

Histo-Line Laboratories

10/06/2023

INDEX

1.Features	3
1.1 Information and specified use.....	3
1.2 Safety and operator comfort.....	4
1.3 Reliability and precision of the cut	5
1.4 Speed and operator safety	7
2.Technical specifications.....	8
2.1 Electrical connections	8
2.2 Section thickness and feed rate	8
2.3 Control features and functions	8
2.4 Dimensions and weight.....	9
2.5 Instrument equipment	9

1.Features

1.1 Information and specified use

S2

The Gemini Histo-Line microtome is a **semi-automatic rotary microtome** particularly suitable for cutting kerosene sections, for normal use in histology, but also for finer cuts used in research. The microtome if necessary can also be used in manual modes.

The microtome is covered with a special solvent-resistant plastic body with no edges, thus minimizing the risk of danger to operators.



The Gemini semiautomatic microtome allows sections to be cut precisely, safely and easily.

It has outstanding ergonomic features for comfortable operation. Gemini is perfectly suited to your cutting habits.

When working with delicate specimens (e.g., brain samples), the microtome will provide reproducible, thin, serial sections of the highest quality. The unique side sample knob and design make sample advancement and retraction more precise. The double handwheel lock system, blade ejector and guard enhance safety.

Digitization of information: Gemini has a 5-inch touch screen that simplifies operation for the operator, changing all parameters with a touch. It has an elegant and simple operating interface that is easy to learn and understand.


1.2 Safety and operator comfort

The flywheel is equipped with ergonomic handle for ease of cutting.

The upper space of the microtome is equipped **with a large area of work surface**, covered with a sheath slide that allows the operator to support the accessories needed for cutting (the package of disposable blades, the brush collect slice, etc.)

The flywheel, which has low friction for good cutting comfort, is equipped with an easily operable dual rotation flywheel lock system through which the sample arm can be stopped in any position. LIBRA is equipped with **two mechanical safety stops**:

- placed on the insertion cutting flywheel.
- the other placed under the toggle flywheel.

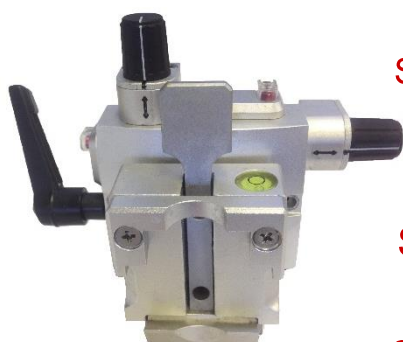

S16 The steel clamp is equipped with a device to "clamp" for quick change of the sample and sample holder fixed tilt.

S1 The clamp is suitable for all types of **standard biocassettes**, and is interchangeable by the operator with the clamp to the biocassettes ring (optional).

S14 The clamp allows easy and precise orientation of specimens with the help of two red **"indicator lights,"** which if clearly visible

indicate that the orientation plane is in the absolute zero position.

On the right side of the clamp is a "bubble" for precise orientation on the three axes of the sample.

Replacing the standard specimen holder clamp is possible through an Allen wrench (included) to use the clamp for super mega cassettes, with a convenient quick coupling.

1.3 Reliability and precision of the cut



The "dual profile" blade holder, suitable for **any type of low- or high-profile disposable blades**, is equipped with the swing-action finger guard. **S10** Through the spring-loaded pin, it is possible to extract or insert the blade **without touching it**.

The blade holder **has three-position** (left, right, center) side-shifting capability to use the entire blade wire without removing the blade. This clearance range corresponds to the length of a standard biocassette so that the entire blade of the microtome can be used without loosening its lock.

The **internal microprocessor** electronically controls the stepper motor so that the cutting and roughing thickness are extremely accurate.

Automatic sample feed is under high-precision electronic control. Through the knobs on the **external control panel**, the thickness and the two speeds of fast and slow feed can be set.

The microtome has the completely **silent sample retraction system suitable** for preserving the life of the blade and the sample, the latter during the return movement is moved to the upper initial position away from the blade. Automatic sample retraction can be excluded by the operator.

Quick and silent return of the sample to the starting position in the roughing (Trimming) stage.

High-precision cutting is suitable for any type of section included in kerosene, resin, routine and laboratory research tasks between 0.5 μm and 100 μm by **stepper motors** and backlash-free, low-maintenance roller bearings.

The Aries microtome has the ability to adjust cutting and roughing thicknesses through the display located on the right side of the microtome or on the external control panel.

The microtome is also designed for extremely fine cutting runs in increments of 0.25 μm to 1 μm



The Aries microtome has a **Touch Screen control display** integrated into the instrument for clear display of information. With this system, it is possible to manage all the functions of the microtome using only the display, which placed on the left side of the operator, gives the indication of the most important working parameters, including:

S16

- cutting thickness
- roughing thickness
- switching of "cutting/roughing" operating modes
- partial and total section counting
- single cut mode
- continuous cutting mode

Both control panels are easy to operate even when wearing protective gloves and can be easily cleaned.

The structure of the microtome is extremely stable, and it is equipped with a powerful motor even for performing particularly hard cuts.

Possibility of using either regrindable blades (by changing the blade holder, an optional accessory) or any type of disposable blades available on the market.

S15

Equipped with an **antistatic tray** for collecting waste sections surrounding the front working area, including the part below the blade holder, which is large, ergonomic for the operator and removable from the tool for easy cleaning.

It is held by two magnets placed in front of the base of the microtome



1.4 Speed and operator safety


The **Aries** and **Gemini** microtomes allow samples to be stored in their sectioning/roughing position and recalled when needed by a sample storage and recognition system. This innovative system can increase operator productivity by automating and therefore reducing the time required to approach each sample to the blade.

To allow rapid exchange of samples (blocks), the instrument is equipped with a "**Memory**" function in which rapid return to the reference cutting position is allowed.

S16

The Memory position allows a custom return position of the holder to be set simply by pressing the button, to rough out samples faster.

The operator can set **two programmable memory positions** by simultaneously pressing two keys located on the display, this allows fast sample repositioning.

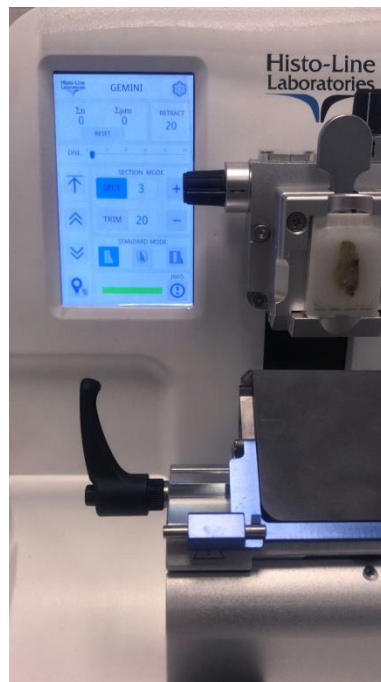

2. Technical specifications

2.1 Electrical connections

Certification: CE, IVD

S3 Power supply: 220/230V~ 50/60Hz

Maximum power consumption: 100W

Fuses: F1A 250V

Safety Class: I

S7 2.2 Section thickness and feed rate

For the "**Cut**" setting, the adjustment range is 0.25-100 µm.

The settable values and cutting increments are:

- 0.25 to 2.5 µm in increments of 0.25 µm
- 2.5 to 5 µm with increments of 0.5 µm
- 5.0 to 10 µm with increments of 1.0 µm
- 10 to 30 µm in increments of 2.0 µm
- 30 to 60 µm in increments of 5.0 µm
- 60 to 100 µm in increments of 10.0 µm

For the "**Trimming**" setting, the adjustment range is 1-600 µm.

The values that can be set and the roughing increments are:

- 1.0 to 10 µm in increments of 1.0 µm
- 10 to 20 µm with increments of 2.0 µm
- 20 to 50 µm in increments of 5.0 µm
- 50 to 150 µm in increments of 10.0 µm
- 150 to 600 µm in increments of 50.0 µm

Retraction specifications:

- S13**
- In manual mode: 5 to 100 µm in increments of 5 µm
 - In automatic mode: speed-dependent, can be deactivated

2.3 Control features and functions

Specimen feed: 30±1 mm, step motor feed

Vertical stroke: 72 mm (suitable for macro sections)

S9 Operating cutting modes: 4 automatic, manual, continuous, single section.

Max. sample size per standard clamp (W x H x D): 55 X 50 X 30 mm

Max. sample block size per clamp for Supermega cassette (H x W x D): 68 x 48 x 15 mm

Sample orientation: both X and Y axes: 0 to +/- 8°, with zero point.

Z axis: 90° rotation.

- S11** Electric rapid feed: 300 $\mu\text{m/s}$ and 800 $\mu\text{m/s}$
Cutting speed: 0.5 to 420 mm/s $\pm 10\%$.
Blade holder displacement adjustment: Forward / Backward: ± 65 mm
Left / Right: ± 40 mm

2.4 Dimensions and weight

Width: 475 mm
Depth (including tray): 580 mm
Depth (excluding tray): 510 mm
Height: 340 mm
Weight: 40 kg

2.5 Instrument equipment

- S1** 1 collection tray sections
1 holder blades
1 clamp for standard samples
Manual use
- S16**

Code: MRS3500