Instruction for use

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«HSV 1,2 IgG-ИФА»

**A SOLID-PHASE ENZYME IMMUNOASSAY
FOR THE QUALITATIVE DETERMINATION
OF IgG ANTIBODIES TO HSV 1, 2
IN HUMAN SERUM OR PLASMA**

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
Authorized Representative in EU:
Polmed.de
Steinacker 20, D-73773
Aichwald, Germany
e-mail: info@polmed.de

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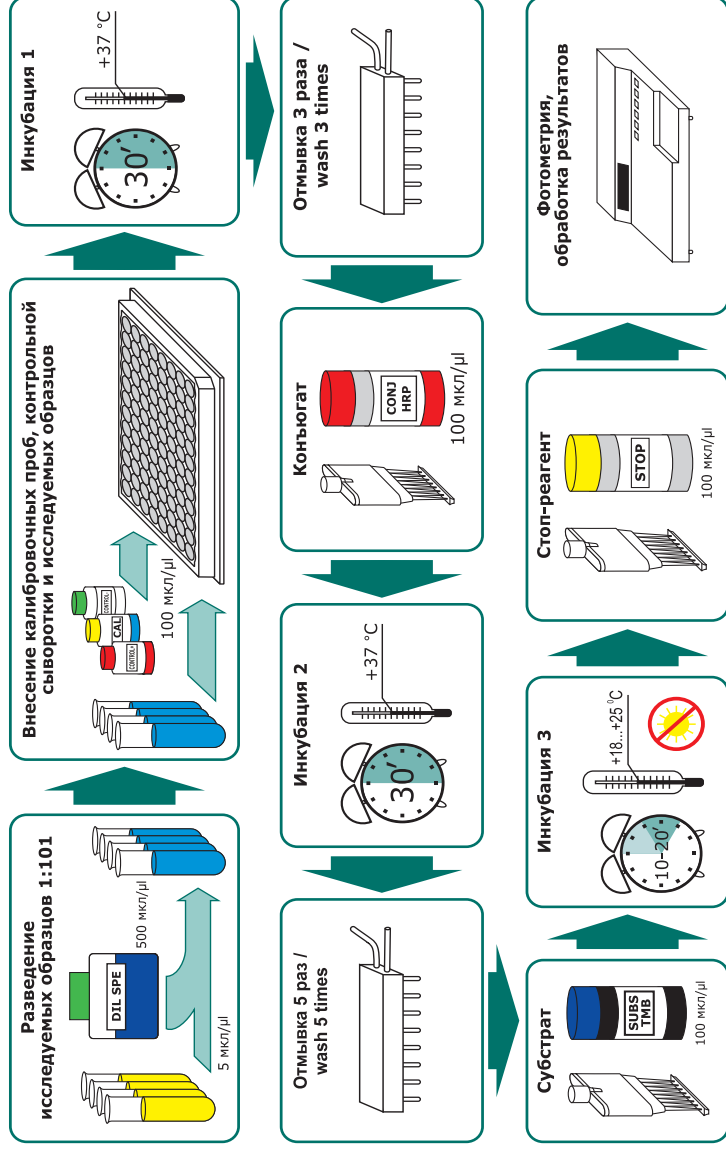CONTENT

1. INTENDED USE	9
2. SUMMARY AND EXPLANATION	9
3. PRINCIPLE OF THE TEST	9
4. WARNINGS AND PRECAUTIONS	10
5. KIT COMPONENTS	11
6. SPECIMEN COLLECTION AND STORAGE	12
7. TEST PROCEDURE	12
8. QUALITY CONTROL	14
9. CALCULATION OF RESULTS	14
10. EXPECTED VALUES	15
11. LITERATURE	15

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Instruction for use

A SOLID-PHASE ENZYME IMMUNOASSAY FOR THE QUALITATIVE DETERMINATION OF IgG ANTIBODIES TO HSV 1,2 IN HUMAN SERUM OR PLASMA

1. INTENDED USE

A solid-phase enzyme immunoassay for the qualitative determination of IgG antibodies to HSV 1,2 in blood serum or plasma.

This kit is designed for measurement of IgG antibodies to HSV 1,2 – in blood serum or plasma. For possibility of use with other sample types, please, refer to Application Notes (on request). The kit contains reagents sufficient for 96 determinations and allows to analyze 45 unknown samples in duplicates.

2. SUMMARY AND EXPLANATION

Herpes simplex virus (HSV) is one of the most common pathogens in humans. HSV is transmitted by all secretions of infected body, especially via saliva, semen and cervical fluid. Acute HSV infections appear as a vesicular rash of labial or genital area. In immunocompromised hosts, HSV may cause life-threatening sequela in central nervous system. HSV is rarely fully eradicated after acute infections and persists in human organism lifelong, showing the periodic reactivation. In case of acute infection or reactivation during pregnancy, HSV may cross the placental barrier and cause severe fetal malformations.

Specific IgG-antibodies to HSV are not protective; their titer usually raise in response to the reactivation of virus and therefore may be used to monitor the actual status of HSV activity. In pregnant women, the absence of specific IgG-antibodies (seronegativity) requires tight restrictions of the lifestyle during the pregnancy minimizing contacts to seropositive humans. The seronegative individuals should not receive the blood transfusions and organ transplants from seropositive donors.

There are two very similar serotypes of HSV – HSV I and HSV II, showing different distribution in affected human tissues and organs. The present test system does not detect the differences between these two serotypes.

3. PRINCIPLE OF THE TEST

This test is based on indirect enzyme immunoassay principle. Tested specimen is placed into the microwells coated by the antigen. Antibodies from the specimen bind coated antigen on the microwell surface. Unbound material is removed by washing procedure. Second antibodies directed towards species specific Ig, labelled with peroxidase enzyme, are then added into the microwells. After subsequent washing procedure, the remaining enzymatic activity bound to the microwell surface is detected and quantified by addition of chromogen-substrate mixture, stop solution and photometry at 450 nm. Optical density in the microwell is directly related to the quantity of the measured analyte in the specimen.

4. WARNINGS AND PRECAUTIONS

4.1. For professional use only.

4.2. This kit is intended for in vitro diagnostic use only.

4.3. INFECTION HAZARD: There is no available test methods that can absolutely assure that Hepatitis B and C viruses, HIV-1/2, or other infectious agents are not present in the reagents of this kit. All human products, including patient samples, should be considered potentially infectious. Handling and disposal should be in accordance with the procedures defined by an appropriate national biohazard safety guidelines or regulations.

4.4. Avoid contact with stop solution containing 5.0% H_2SO_4 . It may cause skin irritation and burns.

4.5. Wear disposable latex gloves when handling specimens and reagents. Microbial contamination of reagents may give false results.

4.6. Do not use the kit beyond the expiration date.

4.7. All indicated volumes have to be performed according to the protocol. Optimal test results are only obtained when using calibrated pipettes and microplate readers.

4.8. Do not smoke, eat, drink or apply cosmetics in areas where specimens or kit reagents are handled.

4.9. Chemicals and prepared or used reagents have to be treated as hazardous waste according to the national biohazard safety guidelines or regulations.

4.10. Do not mix reagents from different lots.

4.11. Replace caps on reagents immediately. Do not swap caps.

4.12. Do not pipette reagents by mouth.

4.13. Specimens must not contain any AZIDE compounds – they inhibit activity of peroxidase.

4.14. Material Safety Data Sheet for this product is available upon request directly from XEMA Co., Ltd.

4.15. The Material Safety Data Sheet fit the requirements of EU Guideline 91/155 EC.

5.1. Contents of the Kit

5. KIT COMPONENTS

Symbol	Description	Qty	Units	Colour code	Stability of opened/diluted components
1 SORB MTP	HSV 1,2 IgG EIA strips, 8x12 wells	1	pcs		until exp.date
2 CAL	Calibrator set, 1.5 ml	1	pcs	blue	2 months
3 CONTROL – CONTROL +	Control sera (1.5 ml)	2	pcs	colourless and red	2 months
4 CONJ HRP	Conjugate, 14 ml	1	pcs	red	until exp.date
5 DIL SPE	EIA buffer 50ml	1	pcs	blue	until exp.date
6 SUBS TMB	Substrate solution, 14 ml	1	pcs	colourless	until exp.date
7 BUF WASH 26X	Washing solution concentrate 26X, 22 ml	1	pcs	colourless	Concentrate – until exp.date Diluted washing solution – 45 days at 2-8 °C or 15 days at RT
8 STOP	Stop solution, 14 ml	1	pcs	colourless	until exp.date
9 N003	Plate sealing tape	2	pcs		N/A
10 K104I	Instruction HSV 1,2 IgG EIA	1	pcs		N/A
11 K104Q	QC data sheet HSV 1,2IgG EIA	1	pcs		N/A

5.2. Equipment and material required but not provided

- Distilled or deionized water;
- Automatic or semiautomatic multichannel micropipettes, 100–250 µl, is useful but not essential;
- Calibrated micropipettes with variable volume, range volume 5–250 µl;
- Calibrated microplate photometer with 450 nm wavelength and OD measuring range 0–3.0;
- Dry thermostat for +37 °C ±0.1 °C.

5.3. Storage and stability of the Kit

Store the whole kit at +2...+8 °C upon receipt until the expiration date.

After opening the pouch keep unused microtiter wells **TIGHTLY SEALED BY ADHESIVE TAPE (INCLUDED)** to minimize exposure to moisture.

6. SPECIMEN COLLECTION AND STORAGE

This kit is intended for use with serum or plasma (ACD- or heparinized). Grossly hemolytic, lipemic, or turbid samples should be avoided.

Specimens may be stored for up to 48 hours at +2...+8 °C before testing. For a longer storage, the specimens should be frozen at -20 °C or lower. Repeated freezing/thawing should be avoided.

7. TEST PROCEDURE**7.1. Reagent Preparation**

- All reagents (including unsealed microstrips) should be allowed to reach room temperature (+18...+25 °C) before use.
- All reagents should be mixed by gentle inversion or vortexing prior to use. Avoid foam formation.
- It is recommended to spin down shortly the tubes with calibrators on low speed centrifuge.
- Prepare washing solution from the concentrate BUF WASH 26X by 26 dilutions in distilled water.

7.2. Procedural Note:

It is recommended that pipetting of all calibrators and samples should be completed within 3 minutes.

7.3. Assay flowchart

See the example of calibration graphic in Quality Control data sheet.

7.4. Assay procedure

1	Put the desired number of microstrips into the frame; allocate 6 wells for the calibrators CAL and control samples CONTROL - , CONTROL + and two wells for each unknown sample. DO NOT REMOVE ADHESIVE SEALING TAPE FROM UNUSED STRIPS.
2	Dilute all samples using buffer DIL SPE (EIA buffer) 101 fold. (5 µl of sample + 500 µl of diluent). Do not dilute control sample and calibrators.
3	Pipet 100 µl of calibrators, control samples CONTROL - , CONTROL + and unknown samples into the wells. Cover the wells by plate adhesive tape (included into the kit).
4	Incubate 30 minutes at +37 °C.
5	Prepare washing solution by 26X dilution of washing solution concentrate BUF WASH 26X by distilled water; Minimal quantity of washing solution should be 250 µl per well. Wash strips 3 times.
6	Dispense 100 µl of CONJ HRP into the wells. Cover the wells by plate adhesive tape.
7	Incubate 30 minutes at +37 °C.
8	Wash the strips 5 times.
9	Dispense 100 µl of SUBS TMB into the wells.
10	Incubate 10-20 minutes at +18...+25 °C.
11	Dispense 100 µl of STOP into the wells.
12	Measure OD (optical density) at 450 nm.
13	Set photometer blank on air.
14	Apply coefficient method for data reduction. 1. Calculate a cutoff value by multiplying mean OD of Calibrator by the Q value (see Lot QC insert). $\text{Cutoff} = \text{OD}(\text{Calibrator}) * Q \text{ value}$ 2. Divide mean OD of each sample by Cutoff. $K = \text{OD}(\text{sample}) / \text{Cutoff}$ 3. If the K value is greater than 1.1, the result is POSITIVE. If the K value is less than 0.9, the result is NEGATIVE. If the K value is between 0.9 and 1.1, the result is EQUIVOCAL.

7.5. Sample processing

Material type	Notes on material collection, storage and handling	Sample dilution example	Sample into the well, µl	Calculation factor
blood serum or plasma	grossly hemolytic, lipemic, or turbid samples should be avoided and should be treated by centrifugation before testing	5 µl of sample + 500 µl of diluent	100	1

8. QUALITY CONTROL

It is recommended to use control samples according to state and federal regulations. The use of control samples is advised to assure the day to day validity of results.

The test must be performed exactly as per the manufacturer's instructions for use. Moreover the user must strictly adhere to the rules of *GLP* (Good Laboratory Practice) or other applicable federal, state, and local standards and/or laws. This is especially relevant for the use of control reagents. It is important to always include, within the test procedure, a sufficient number of controls for validating the accuracy and precision of the test.

The test results are valid only if all controls are within the specified ranges and if all other test parameters are also within the given assay specifications.

9. CALCULATION OF RESULTS

Some laboratories, based on their population studies, set up a second cutoff, which stands between anamnestic ('normal') IgG antibody level and 'high' IgG antibody level characteristic of reactivation or late period of primary infection. Recommended values for this second cut-off for two age groups (8 months – 3 year, > 3 years) are presented in the table below.

10. EXPECTED VALUES

Therapeutical consequences should not be based on results of IVD methods alone – all available clinical and laboratory findings should be used by a physician to elaborate therapeutically measures. Each laboratory should establish its own normal range for HSV 1,2 IgG. Based on data obtained by XEMA, the following normal range is recommended (see below).

Sex, age	Units, K	
	Lower limit	Upper limit
Seronegative	<0.1	0.9
Seropositive > 3 years	1.1	4.9
newborn*	<0.1	1.5
under 8 months*	<0.1	3.5
8 months – 3 years	<0.1	4.8

*antibodies of maternal origin

11. LITERATURE

1. Isakov V. A., Borisova V. V., Isakov D. V. Herpetic infection: pathogenesis and laboratory diagnostics. Manual for physicians. St.-Petersburg, Lanj publishing, 1999, p. 192.
2. Cowan F. M., Johnson A. M., Ashley R. et al. Relation ship between antibodies to herpes simplex virus (HSV) and symptoms of HSV infection // J. Infect. Dis. 1996- v.174 p.470-475
3. Munday P. E., Vuddamalay J., Slomka M. J. et al. Role of type specific herpes simplex virus serilogy in the diagnosis and management of genital herpes // Sex Transm. Inf 1998 v.74 p.175-178
4. Roest R. W., van der Meijden W. I., van Dijk G. et al. Prevalence and association between herpes simplex viruses types 1- and 2-specific antibodies in attendees at a sexually transmitted disease clinic. // Int Epidemiol. 2001 v.30 p.580-588
5. Whitley R/I/ Herpes simplex virusinfections. In: Infection Diseases of the Fetus and Newborn Infants (Remington JS, Klein JO eds.) Philadelphia: WB Saunders, 1990 - p.282-305

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Russian Diagnostic
Manufacturers Association

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RUSSIAN ASSOCIATION
OF MEDICAL LABORATORY
DIAGNOSTICS

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