Instruction for use

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**A SOLID-PHASE ENZYME IMMUNOASSAY
FOR THE QUANTITATIVE DETERMINATION
OF FSH IN HUMAN SERUM OR PLASMA**

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
e-mail: redkin@xema-medica.com

internet: www.xema-medica.com

Authorized Representative in EU:
Polmed.de
Steinacker 20, D-73773
Aichwald, Germany
e-mail: info@polmed.de

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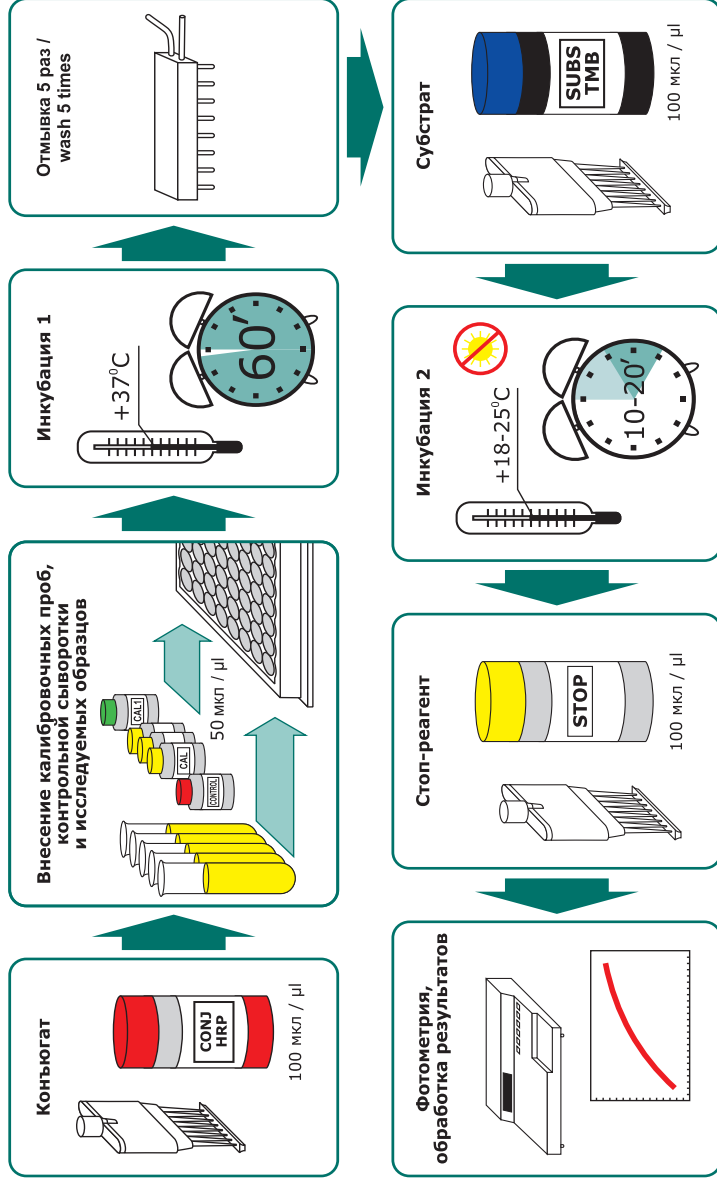CONTENT

1. INTENDED USE	10
2. SUMMARY AND EXPLANATION	10
3. PRINCIPLE OF THE TEST	11
4. WARNINGS AND PRECAUTIONS	11
5. KIT COMPONENTS	12
6. SPECIMEN COLLECTION AND STORAGE	13
7. TEST PROCEDURE	13
8. QUALITY CONTROL	14
9. CALCULATION OF RESULTS	15
10. EXPECTED VALUES	15
11. PERFORMANCE CHARACTERISTICS	16
12. LITERATURE	16

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Instruction for use

A SOLID-PHASE ENZYME IMMUNOASSAY FOR THE QUANTITATIVE DETERMINATION OF FSH IN HUMAN SERUM OR PLASMA

1. INTENDED USE

A solid-phase enzyme immunoassay for the quantitative determination of FSH in blood serum or plasma.

This kit is designed for measurement of FSH in blood serum or plasma. For possibility of use with other sample types, please, refer to Application Notes (on request). The kit contains reagents sufficient for 96 determinations and allows to analyze 42 unknown samples in duplicates.

2. SUMMARY AND EXPLANATION

Follicle stimulating hormone (FSH) is a glycoprotein with molecular weight 28 kDa secreted by basophil cells in hypophysis. Gonadotropin releasing hormone (GnRH) produced by the hypothalamus controls the release of FSH from anterior pituitary. Follicle-stimulating hormone (FSH) and Luteinizing hormone (LH) are intimately involved in the control of the growth and reproductive activities of the gonadal tissues, which synthesize and secrete male and female sex hormones. The levels of circulating FSH and LH are controlled by these sex hormones through a negative feedback. Like LH, TSH and HCG, FSH consists of two subunits – alpha and beta, its biological and immunological properties being dependent on the hormone-specific beta subunit.

In females, FSH stimulates the growth and maturation of ovarian follicles. At the beginning of normal menstrual cycle FSH level is higher than at the final stage of follicular phase. Peak FSH levels are seen in the middle of the cycle concomitantly with LH peak levels. Increased estradiol and progesterone production during luteal phase leads to decreased FSH blood concentrations by negative feedback mechanism. The same mechanism leads to elevation of FSH levels at the end of the cycle due to decreased estrogen and progesterone concentrations, and the new cycle is initiated.

In men, FSH regulates growth of seminiferous tubules and maintenance of spermatogenesis. However, androgens, unlike estrogen, do not lower FSH level, therefore demonstrating a feedback relationship only with serum LH. High levels of FSH in women are seen in menopause, preliminary ovarian failure, agenesis of ovaries; in men elevated FSH levels may be found in primary testicular failure, dysgenesis of seminiferous tubules, delayed sexual maturation, and Klinefelter syndrome. Elevated concentrations are also found in cases of starvation, renal failure, hyperthyroidism, cirrhosis and after intake of clomifen, L-DOPA.

Decreased FSH levels are found in hypopituitarism and after intake of oral contraceptives, phenothiazine, estrogens.

3. PRINCIPLE OF THE TEST

This test is based on two-site sandwich enzyme immunoassay principle. Tested specimen is placed into the microwells coated by specific murine monoclonal to β chain human FSH-antibodies. Antigen from the specimen is captured by the antibodies coated onto the microwell surface. Second antibodies - murine monoclonal to α chain human LH/FSH/HCG, labelled with peroxidase enzyme, are then added into the microwells. After washing procedure, the remaining enzymatic activity bound to the microwell surface is detected and quantified by addition of chromogen-substrate mixture, stop solution and photometry at 450 nm. Optical density in the microwell is directly related to the quantity of the measured analyte in the specimen.

4. WARNINGS AND PRECAUTIONS

4.1. For professional use only.

4.2. This kit is intended for in vitro diagnostic use only.

4.3. INFECTION HAZARD: There is no available test methods that can absolutely assure that Hepatitis B and C viruses, HIV-1/2, or other infectious agents are not present in the reagents of this kit. All human products, including patient samples, should be considered potentially infectious. Handling and disposal should be in accordance with the procedures defined by an appropriate national biohazard safety guidelines or regulations.

4.4. Avoid contact with stop solution containing 5.0% H_2SO_4 . It may cause skin irritation and burns.

4.5. Wear disposable latex gloves when handling specimens and reagents. Microbial contamination of reagents may give false results.

4.6. Do not use the kit beyond the expiration date.

4.7. All indicated volumes have to be performed according to the protocol. Optimal test results are only obtained when using calibrated pipettes and microplate readers.

4.8. Do not smoke, eat, drink or apply cosmetics in areas where specimens or kit reagents are handled.

4.9. Chemicals and prepared or used reagents have to be treated as hazardous waste according to the national biohazard safety guidelines or regulations.

4.10. Do not mix reagents from different lots.

4.11. Replace caps on reagents immediately. Do not swap caps.

4.12. Do not pipette reagents by mouth.

4.13. Specimens must not contain any AZIDE compounds – they inhibit activity of peroxidase.

4.14. Material Safety Data Sheet for this product is available upon request directly from XEMA Co., Ltd.

4.15. The Material Safety Data Sheet fit the requirements of EU Guideline 91/155 EC.

5.1. Contents of the Kit

5. KIT COMPONENTS

Symbol	Description	Qty	Units	Colour code	Stability of opened/diluted components
1 SORB MTP	FSH EIA strips, 8x12 wells	1	pcs		until exp.date
2 CAL 1 - 5	polystyrene microwells coated with murine monoclonal to β chain human FSH	5	pcs	green (C1 - colourless)	2 months
3 CONTROL	human FSH diluted in a preselected animal serum preservative - 0.01% Bronidox L, 0.01% 2-Methyl-4-isothiazolin-3-one-hydrochloride; also contains green dye	1	pcs	colourless	2 months
4 CONJ HRP	dilution of preselected human serum, with high content of FSH with casein solution; preservative - 0.1% phenol, colourless	1	pcs	green	until exp.date
5 SUBS TMB	aqueous solution of murine monoclonal to a chain human LH/FSH/HCG coupled with horseradish peroxidase diluted on phosphate buffered solution with casein from bovine milk and detergent (Tween-20), contains 0.1% phenol as preservative and green dye	1	pcs	colourless	until exp.date
6 BUF WASH 26X	ready-to-use single-component tetramethylbenzidine (TMB) solution.	1	pcs	colourless	Concentrate - until exp.date Diluted washing solution - 45 days at 2-8 °C or 15 days at RT
7 STOP	aqueous solution of sodium chloride and detergent (Tween 20), contains proClin300 as a preservative	1	pcs	colourless	until exp.date
8 N003	5.0% vol/vol solution of sulphuric acid	2	pcs		N/A
9 K2031	Stop solution, 14 ml Plate sealing tape Instruction FSH EIA	1	pcs		N/A
10 K203Q	QC data sheet FSH EIA	1	pcs		N/A

5.2. Equipment and material required but not provided

- Distilled or deionized water;
- Automatic or semiautomatic multichannel micropipettes, 100–250 µl, is useful but not essential;
- Calibrated micropipettes with variable volume, range volume 25–250 µl;
- Dry thermostat for +37 °C ±0.1 °C;
- Calibrated microplate photometer with 450 nm wavelength and OD measuring range 0–3.0

5.3. Storage and stability of the Kit

Store the whole kit at +2...+8 °C upon receipt until the expiration date.

After opening the pouch keep unused microtiter wells TIGHTLY SEALED BY ADHESIVE TAPE (INCLUDED) to minimize exposure to moisture.

6. SPECIMEN COLLECTION AND STORAGE

This kit is intended for use with serum or plasma (ACD- or heparinized). Grossly hemolytic, lipemic, or turbid samples should be avoided.

Specimens may be stored for up to 48 hours at +2...+8 °C before testing. For a longer storage, the specimens should be frozen at -20 °C or lower. Repeated freezing/thawing should be avoided.

7. TEST PROCEDURE**7.1. Reagent Preparation**

- All reagents (including unsealed microstrips) should be allowed to reach room temperature (+18...+25 °C) before use.
- All reagents should be mixed by gentle inversion or vortexing prior to use. Avoid foam formation.
- It is recommended to spin down shortly the tubes with calibrators on low speed centrifuge.
- Prepare washing solution from the concentrate BUF WASH 26X by 26 dilutions in distilled water.

7.2. Procedural Note:

It is recommended that pipetting of all calibrators and samples should be completed within 3 minutes.

7.3. Assay flowchart

See the example of calibration graphic in Quality Control data sheet.

7.4. Assay procedure

1	Put the desired number of microstrips into the frame; allocate 12 wells for the calibrators CAL 1–5 and control samples CONTROL and two wells for each unknown sample. DO NOT REMOVE ADHESIVE SEALING TAPE FROM UNUSED STRIPS.
2	If suggested analyte concentration in the sample exceeds the highest calibrator, additionally dilute this sample accordingly, using (zero calibrator). Use of other buffers or reagents for sample dilution may lead to incorrect measurement.
3	Dispense 100 µl of CONJ HRP into the wells.
4	Pipet 50 µl of calibrators CAL 1–5, control samples CONTROL and unknown samples into the wells. Cover the wells by plate adhesive tape (included into the kit).
5	Incubate 60 minutes at 37 °C.
6	Prepare washing solution by 26X dilution of washing solution concentrate (BUF WASH 26X) by distilled water. Wash the strips 5 times.
7	Dispense 100 µl of SUBS TMB into the wells
8	Incubate 10–20 minutes at +18...+25 °C
9	Dispense 100 µl of STOP into the wells.
10	Measure OD (optical density) at 450 nm.
11	Set photometer blank on first calibrator
12	Apply point-by-point method for data reduction.

7.5. Handling notes

Calibrators and control sample(s) – only one freezing/thawing cycle is allowed

8. QUALITY CONTROL

It is recommended to use control samples according to state and federal regulations. The use of control samples is advised to assure the day to day validity of results.

The test must be performed exactly as per the manufacturer's instructions for use. Moreover the user must strictly adhere to the rules of GLP (Good Laboratory Practice) or other applicable federal, state, and local standards and/or laws. This is especially relevant for the use of control reagents. It is important to always include, within the test procedure, a sufficient number of controls for validating the accuracy and precision of the test.

The test results are valid only if all controls are within the specified ranges and if all other test parameters are also within the given assay specifications.

9. CALCULATION OF RESULTS


9.1. Calculate the mean absorbance values (OD450) for each pair of calibrators and samples.

9.2. Plot a calibration curve on graph paper: OD versus FSH concentration.

9.3. Determine the corresponding concentration of FSH in unknown samples from the calibration curve. Manual or computerized data reduction is applicable on this stage. Point-by-point or linear data reduction is recommended due to non-linear shape of curve.

9.4. Below is presented a typical example of a standard curve with the XEMA Co. Not for calculations!

Calibrators	Value	Absorbance Units (450 nm)
CAL 1	0 IU/l	0.05
CAL 2	5 IU/l	0.27
CAL 3	25 IU/l	0.87
CAL 4	50 IU/l	1.40
CAL 5	100 IU/l	2.41

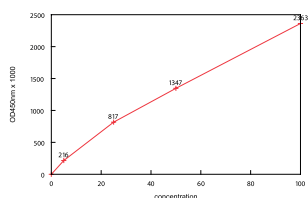10. EXPECTED VALUES

Therapeutical consequences should not be based on results of IVD methods alone – all available clinical and laboratory findings should be used by a physician to elaborate therapeutically measures. Each laboratory should establish its own normal range for FSH. Based on data obtained by XEMA, the following normal range is recommended (see below). NOTE: the patients that have received murine monoclonal antibodies for radioimaging or immunotherapy develop high titered anti-mouse antibodies (HAMA). The presence of these antibodies may cause false results in the present assay. Sera from HAMA positive patients should be treated with depleting adsorbents before assaying.

Sex, age	Units, IU/l	
	Lower limit	Upper limit
Children under 11 yrs	-	4.0
Males	0.8	25
Females		
Menstrual cycle:		
follicular phase	3.0	12
luteinic phase	2.0	12
ovulation	6.0	25
post menopausal	10.0	150

11. PERFORMANCE CHARACTERISTIC**11.1.** Analytical specificity / Cross reactivity

Analyte	Cross-reactivity, % wt/wt
TSH	<0.1
LH	<0.1
HCG	<0.1

11.2. Analytical sensitivity

Sensitivity of the assay was assessed as being 0.15 IU/l.

11.3. Linearity


Linearity was checked by assaying dilution series of 5 samples with different FSH concentrations. Linearity percentages obtained ranged within 90 to 110%.

11.4. Recovery

Recovery was estimated by assaying 5 mixed samples with known FSH concentrations. The recovery percentages ranged from 90 to 110%.

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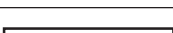Russian Diagnostic
Manufacturers Association

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Russian Association
of Medical Laboratory
Diagnosticians

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