

QuantStudio™ Design and Analysis Desktop Software

USER GUIDE

Getting started with design and analysis of experiments in the desktop software v1.6.

for use with:

QuantStudio™ 1 Real-Time PCR System

QuantStudio™ 3 Real-Time PCR System

QuantStudio™ 5 Real-Time PCR System

Publication Number MAN0010408

Revision E.0



For Research Use Only. Not for use in diagnostic procedures.

ThermoFisher
S C I E N T I F I C



Revision history: MAN0010408 E.0 (English)

Revision	Date	Description
E.0	30 January 2024	v1.6 • Network and password security requirements were added. • A chapter was added to describe the use of the software with the Security, Auditing, and E-signature Administrator Console. • Additional instructions were provided for exporting run data and results ("Export experiments or results" on page 42). • The preferences were updated to include selecting the data destination for different file types ("Set global desktop software preferences" on page 15). This replaces selecting a single default folder or selecting the last location where a file was saved. • The instructions to print a report were updated to note that the grouping in the well table is retained in the report ("Print a report of run data and results" on page 44). • Instructions to save as were updated to include a new dialog box. • The required permissions were included in the applicable sections when security is enabled in the software.
D.0	22 May 2023	The references to the specific user guides for each real-time PCR instrument were removed.
C.0	19 October 2018	Updated to add QuantStudio™ 1 Real-Time PCR System information.
B.0	17 December 2015	Updates include: • Combination of define and assign functions into a single Plate tab • Display of VeriFlex™ Zones on plate layout • Real-time data monitoring in the Run tab • Security, Audit, and E-Signature (SAE) features • Implementation of locked workflow • Selection of an instrument before starting a run • Ability to select multiple targets in the results view • Various minor changes to the user interface
A.0	17 April 2015	New document.

The information in this guide is subject to change without notice.

DISCLAIMER: TO THE EXTENT ALLOWED BY LAW, THERMO FISHER SCIENTIFIC INC. AND/OR ITS AFFILIATE(S) WILL NOT BE LIABLE FOR SPECIAL, INCIDENTAL, INDIRECT, PUNITIVE, MULTIPLE, OR CONSEQUENTIAL DAMAGES IN CONNECTION WITH OR ARISING FROM THIS DOCUMENT, INCLUDING YOUR USE OF IT.

NOTICE TO PURCHASER: DISCLAIMER OF LICENSE: Purchase of this software product alone does not imply any license under any process, instrument or other apparatus, system, composition, reagent or kit rights under patent claims owned or otherwise controlled by Thermo Fisher Scientific, either expressly, or by estoppel.

TRADEMARKS: All trademarks are the property of Thermo Fisher Scientific and its subsidiaries unless otherwise specified. TaqMan is a trademark of Roche Molecular Systems, Inc., used under permission and license. Microsoft, Excel, and PowerPoint are trademarks of Microsoft Corporation.

©2015-2024 Thermo Fisher Scientific Inc. All rights reserved.

Contents

■	CHAPTER 1	Product information	9
		Product description	9
		Network and password security requirements	9
		Network configuration and security	9
		Password security	9
		Recommendations for passwords	10
		Features in the QuantStudio™ Design and Analysis Desktop Software	10
		Experiment types	11
		Workflow overview	14
■	CHAPTER 2	General procedures to set up and run experiments	15
		Set global desktop software preferences	15
		Workflow: Set up and run an experiment	16
		Set up a template	16
		Create or open a template	17
		Enter template properties	17
		Confirm or edit the run method and optical filter selection	19
		Assign plate and well attributes using the Quick Setup subtab	20
		Save a template file	22
		Prepare reactions	23
		Good laboratory practices for PCR and RT-PCR	23
		Guidelines for handling samples and reagents	23
		Guidelines for setting up the reactions in the plates or tubes	24
		Start and monitor a run	24
		View real-time run information in the desktop software	25
		View real-time run information on the instrument touchscreen	25
		View the Post-Run Summary in the desktop software	26
		Experiment setup using libraries, external files, and templates	27
■	CHAPTER 3	General procedures to review results	28
		About the quantification cycle (C _q)	28
		Overview of the Results tab	29
		Workflow: General procedures to review the run results	31
		Guidelines for viewing and analyzing results	31

Assess results in the Amplification Plot	32
Amplification Plot overview	32
Assess the overall shape of the Amplification Plot curves	32
Confirm or correct threshold settings	33
Confirm or correct baseline settings	34
Omit outliers from analysis	35
Optimize display of negative controls in the Amplification Plot	36
Assess results in the Well Table view	37
Well Table overview	37
Group or sort the Well Table	37
Review the dye signal profile using the Multicomponent Plot	38
Multicomponent Plot overview	38
View and assess the Multicomponent Plot	38
Review the signal profile using the Raw Data Plot	39
Raw Data Plot overview	39
View and assess the Raw Data Plot	40
Review the flags in the QC Summary	40
View calibration results in the desktop software	41
Review ROI/Uniformity calibration results	41
Review Background calibration results	41
Review Dye calibration results	41
Override the calibration data	42
Export experiments or results	42
Options for exporting the plots	42
Export run data results	43
Print the plate layout	44
Send run results to PowerPoint™	44
Print a report of run data and results	44
Export configurations	45
 ■ CHAPTER 4 Set up, run, and review standard curve experiments	46
Standard curve experiments	46
Overview	46
Reaction types	46
Compatible PCR options	47
Set up a standard curve experiment in the software	47
Set up and run the PCR reactions	48
Review results	49
Workflow: Review standard curve experiments	49
Standard Curve Plot overview	50
View and assess the Standard Curve Plot	51
Standard curve settings overview	52

■ CHAPTER 5	Set up, run, and review	
	relative standard curve experiments and comparative C_t experiments	53
	Relative standard curve experiments	53
	Overview	53
	Reaction types	54
	Compatible PCR options	54
	Comparative C_t experiments	55
	Overview	55
	Reaction types	55
	Compatible PCR options	56
	Relative quantitation: relative standard curve vs. comparative C_t	56
	Set up a relative standard curve experiment in the software	57
	Set up a comparative C_t experiment in the software	58
	Set up and run the PCR reactions	59
	Review results	60
	Workflow: Review relative standard curve and comparative C_t experiments	60
	View and assess the Standard Curve Plot	61
	Gene Expression Plot overview	63
	QC Plot overview	64
	Relative quantification settings overview	65
■ CHAPTER 6	Set up, run, and review genotyping experiments	66
	Genotyping experiments	66
	Overview	66
	Reaction types	67
	Compatible PCR options	67
	Set up a genotyping experiment in the software	67
	Set up and run the PCR reactions	68
	Review results	69
	Workflow: Review genotyping experiments	69
	Allelic Discrimination Plot overview	69
	View and assess the Allelic Discrimination Plot	70
	Perform manual calls	71
	Call settings overview (genotyping)	71

■	CHAPTER 7	Set up, run, and review presence/absence experiments	73
	Presence/absence experiments		73
	Overview		73
	Reaction types		73
	Compatible PCR options		74
	Set up the presence/absence experiment in the software		74
	Set up and run the PCR reactions		75
	Review results		76
	Workflow: Review presence/absence experiments		76
	Presence/Absence Plot overview		76
	View and assess the Presence/Absence Plot		77
	Call settings overview (presence / absence)		78
■	CHAPTER 8	Set up, run, and review melt curve experiments	80
	Melt curve experiments		80
	Overview		80
	Reaction types		80
	Set up a melt curve experiment in the software		81
	Set up melt curve reactions		82
	Review results		82
	Workflow: Review melt curve experiments		82
	Melt Curve Plot overview		82
	View and assess the Melt Curve Plot		83
	Melt curve settings overview		83
■	CHAPTER 9	Use the software with the Security, Auditing, and E-signature (SAE) Administrator Console	85
	Overview of the Security, Auditing, and E-signature Administrator Console		85
	Enable the system security function		86
	Disable the system security function		86
	Overview of accounts		87
	Log in to QuantStudio™ Design and Analysis Desktop Software		87
	Change your password in QuantStudio™ Design and Analysis Desktop Software		87
	Log out of QuantStudio™ Design and Analysis Desktop Software		88
	Default roles and permissions		88

View and manage the audit records	89
View the object audit records	89
Generate an object audit report	90
Export the audit summary	90
Export the audit details	91
View and manage the e-signatures	91
Sign the data	91
View the e-signatures	91
■ APPENDIX A Instrument overview	93
Power on the instrument	93
View run history and delete or transfer files from the instrument	93
Load and unload the plate in the instrument	94
Load and unload a plate in the QuantStudio™ 1 Real-Time PCR Instrument	94
Load and unload a plate in the QuantStudio™ 3 Real-Time PCR Instrument or QuantStudio™ 5 Real-Time PCR Instrument	95
■ APPENDIX B Alternative procedures to set up a template	97
Set up a custom experiment	97
Custom experiments overview	97
Set up a custom experiment in the software	98
Assign samples using a sample definition file	99
About sample definition files	99
Create a sample definition file	99
Import sample information from a sample definition file	100
Assign samples and targets using plate setup files	100
About plate setup files	100
Import plate setup data	100
Assign targets, samples, and biological replicate groups from an XLS file	101
Create new EDT files using existing EDT and EDS files	102
About experiment templates	102
■ APPENDIX C Detailed procedures to create or edit a method	103
Adjust method parameters	103
Set up advanced temperature zones (Auto Delta and VeriFlex™ Zones)	104
Add or adjust a pause step	105
Select optical filters	105
Edit the ramp increment for the melt curve dissociation step	106

■ APPENDIX D Detailed procedures to set up plate / well details and libraries	107
Assign well attributes (Quick Setup subtab)	107
Define and assign well attributes (Advanced Setup subtab)	107
Assign a task to wells	108
Detection tasks for targets and SNP assays	108
Define and set up standard dilutions	109
Assign the standard dilutions manually	110
Define and assign biological replicate groups	110
Sample, target, and SNP assay libraries	111
Libraries overview	111
Permissions for the libraries	111
Apply a filter to search a library	111
Define samples in the Samples table	112
Define targets in the Targets table	113
Define SNP assays in the SNP Assays table	114
■ APPENDIX E Configure analysis settings	115
Guidelines for the analysis settings	115
View and configure the analysis settings	116
C _t settings overview	116
Flag settings overview	117
Advanced settings overview	117
■ APPENDIX F Documentation and support	118
Related documentation for the QuantStudio™ 1 Real-Time PCR System	118
Related documentation for the QuantStudio™ 3 Real-Time PCR System and QuantStudio™ 5 Real-Time PCR System	118
Customer and technical support	118
Limited product warranty	119
Glossary	120
Index	122



Product information

■ Product description	9
■ Network and password security requirements	9
■ Features in the QuantStudio™ Design and Analysis Desktop Software	10
■ Experiment types	11
■ Workflow overview	14

Product description

The QuantStudio™ Design and Analysis Desktop Software allows the user to open, run, and analyze experiments generated with QuantStudio™ 1 Real-Time PCR System, QuantStudio™ 3 Real-Time PCR System, and QuantStudio™ 5 Real-Time PCR System. The software also allows you to set up experiments, send experiments to the instrument, collect data, and analyze the collected data.

The software can be used with the Security, Auditing, and E-signature Administrator Console. The console enables system security by controlling user access, facilitating auditing, and supporting e-signatures. For more information, see Chapter 9, “Use the software with the Security, Auditing, and E-signature (SAE) Administrator Console”.

Network and password security requirements

Network configuration and security

The network configuration and security settings of your laboratory or facility (such as firewalls, anti-virus software, network passwords) are the sole responsibility of your facility administrator, IT, and security personnel. This product does not provide any network or security configuration files, utilities, or instructions.

If external or network drives are connected to the software, it is the responsibility of your IT personnel to ensure that such drives are configured and secured correctly to prevent data corruption or loss. It is the responsibility of your facility administrator, IT, and security personnel to prevent the use of any unsecured ports (such as USB, Ethernet) and ensure that the system security is maintained.

Password security

Thermo Fisher Scientific strongly recommends that you maintain unique passwords for all accounts in use on this product. All passwords should be reset upon first sign in to the product. Change passwords according to your organization's password policy.

It is the sole responsibility of your IT personnel to develop and enforce secure use of passwords.

Recommendations for passwords

Thermo Fisher Scientific recommends enabling a password policy for SAE user accounts with the following minimum number of characters.

- Administrative users—12 characters
- Non-administrative users—8 characters

The use of a password manager is recommended to help create secure passwords.

Features in the QuantStudio™ Design and Analysis Desktop Software

Actions	Unlocked template	Locked template ^[1,2]
Properties tab		
Edit experiment name; enter / scan plate barcode; enter user name	✓	✓
Select instrument / block type, experiment type, chemistry, run mode	✓	—
Enter reagent information (chemistry details)	✓	✓
Method tab		
Edit the thermal protocol, reaction volume, optical filter selection	✓	—
Plate tab		
Quick Setup subtab		
Define plate attributes	✓	—
Define or assign samples	✓	✓
Define or assign targets or SNP assays	✓	—
Advanced Setup subtab		
Define samples	✓	✓
Assign samples	✓	✓
Define targets or SNP assays	✓	—
Assign targets or SNP assays	✓	✓
Run tab		
Start and monitor a run in progress	✓	✓
View time remaining and plots ^[3]	✓	✓

(continued)

Actions	Unlocked template	Locked template ^[1,2]
Results tab		
Review run results (analyzed run data)	✓	✓
Configure analysis settings	✓	—
Export tab		
Select export options for run data and run results (analyzed run data)	✓	✓
Export run data ^[3] ; export template settings	✓	✓

^[1] If you enter the password for a locked template, all listed actions are available.

^[2] Always save a backup *unlocked* version of a template before saving it as a locked template.

^[3] This feature is also available from the instrument touchscreen.

Experiment types



Purpose	Description
Standard curve experiment	
Determines absolute target quantity in samples.	1. The software measures amplification of the target in a standard dilution series and in test samples. 2. The software generates a standard curve using data from the standard dilution series. 3. The software uses the standard curve to interpolate the absolute quantity of target in the test samples.

Purpose	Description
Relative standard curve experiment	
Determines relative target quantity in samples.	1. The software measures amplification of the target of interest and of an endogenous control target in a standard dilution series, in a reference (calibrator) sample, and in test samples. The endogenous control is a target that is expressed equally in all samples; examples of endogenous controls are β-actin, GAPDH, and 18S ribosomal RNA. The software can algorithmically incorporate multiple endogenous control targets in relative quantification calculations. The reference sample is used as the basis for relative quantification results (or 1× sample). For example, in a study of drug effects on gene expression, an untreated control is an appropriate reference sample. 2. The software generates standard curves for the target of interest and the endogenous control using data from the corresponding standard dilution series. 3. The software uses the standard curves to interpolate the quantities of the target of interest and the endogenous control in each sample. The target quantity in each sample is then normalized to the sample's endogenous control quantity. 4. To determine the relative quantity of the target in test samples, the software divides the normalized target quantity in the sample by the normalized target quantity in the reference sample.
Comparative C_t ($\Delta\Delta C_t$) experiment	
Determines relative target quantity in samples.	1. The software measures amplification of the target of interest and of an endogenous control target in a reference (calibrator) sample and in test samples. The endogenous control is a target that is expressed equally in all samples; examples of endogenous controls are β-actin, GAPDH, and 18S ribosomal RNA. The software can algorithmically incorporate multiple endogenous control targets in relative quantification calculations. The reference sample is used as the basis for relative quantification results (or 1× sample). For example, in a study of drug effects on gene expression, an untreated control is an appropriate reference sample. 2. The measurements for the target of interest are normalized to the endogenous control. 3. To determine the relative quantity of the target in test samples, the software compares the normalized ΔC_q (ΔC_t or ΔC_{rt}) for the sample to the normalized ΔC_q (ΔC_t or ΔC_{rt}) for the reference sample.

Purpose	Description
Genotyping experiment	
Detects single nucleotide polymorphism (SNP) variants of a target nucleic acid sequence.	Genotyping experiments use preformulated TaqMan™ SNP Genotyping Assays that include the following components: • Two sequence-specific primers for amplification of sequences containing the SNP of interest • Two allele-specific TaqMan™ probes for Allele 1 and Allele 2 1. The software normalizes the fluorescence of the reporter dyes to the fluorescence of the passive reference dye in each well. 2. The software plots the normalized reporter dye signal of each sample well on an Allelic Discrimination Plot, which contrasts the reporter dye intensities of the allele-specific probes. 3. The software algorithmically clusters the sample data, and assigns a genotype call to the samples of each cluster according to its position on the plot.
Presence/absence experiment	
Determines the presence or absence of a target nucleic acid sequence in a sample.	The software calls the target present or absent based on an algorithmically determined call threshold. (The call threshold is different from the C_t threshold; the C_t threshold is not used to make calls.)
Melt curve experiment	
Determines the melting temperature (T_m) of the amplification products of a PCR that used intercalating dyes.	In the software, melt curve analysis is included in the default run method for any experiment type that uses intercalating dyes. 1. The software plots a melt curve based on the fluorescence of the dye with respect to change in temperature. 2. Using the melt curve, the software calculates the melting temperature (T_m).

Workflow overview

Use the tabs across the top of the screen to navigate the workflow in the QuantStudio™ desktop software.

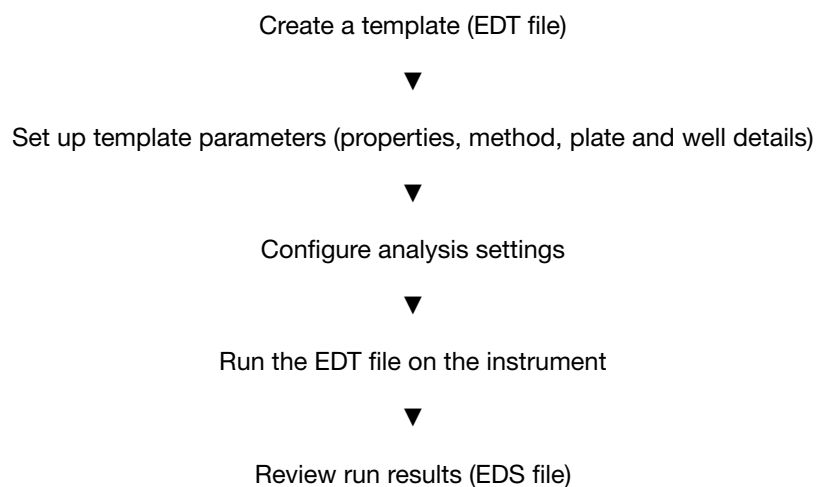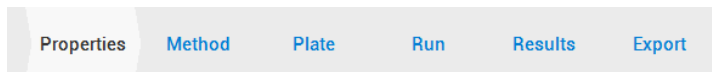

2

General procedures to set up and run experiments

■ Set global desktop software preferences	15
■ Workflow: Set up and run an experiment	16
■ Set up a template	16
■ Prepare reactions	23
■ Start and monitor a run	24
■ Experiment setup using libraries, external files, and templates	27

Set global desktop software preferences

If security is enabled in the software, your account must have the permission of **Edit System Preferences**.

1. In the menu bar, select **Tools ▶ Preferences**.
2. In the **System Preferences** dialog box, set preferences in each tab.

Tab	Options	Action
Defaults	Decimal Places to Show	Enter the number of significant figures in exported results.
	Language [if available]	Choose a language for the software from the dropdown menu.
Experiment	Baseline Start Cycle and Baseline End Cycle	Enter the first and last cycles to be used to calculate the baseline for runs that include amplification.
	Auto Analysis	Select to perform auto analysis at the end of each run.
	Auto Save	Select to save changes at the end of each run.
Print	Disable color when printing the Well Table	Select to disable color printing.
Data Destination	Save As Location, Export Location, and Report Location	Select where to save files and where to export results. Click Browse , then navigate to and select a default location.
Display Format	Date Format, Time Format, and Decimal Point Format	Select the display formats. These formats are also used in the export or import of data.

3. Click **Save**.

Workflow: Set up and run an experiment

Set up a template (page 16)

Create or open a template (page 17)

Enter template properties (page 17)

Confirm or edit the run method and optical filter selection (page 19)

Assign plate and well attributes (page 20)

Save a template file (page 22)



Prepare reactions (page 23)



Start and monitor a run (page 24)

For post-run procedures, see the following sections:

- Chapter 3, “General procedures to review results”
- “Export experiments or results” on page 42

Set up a template

This section describes the general procedures to set up a template in the desktop software. For setup information for a specific experiment type, see the following sections.

- “Set up a standard curve experiment in the software” on page 47
- “Set up a relative standard curve experiment in the software” on page 57
- “Set up a comparative Ct experiment in the software” on page 58
- “Set up a genotyping experiment in the software” on page 67
- “Set up the presence/absence experiment in the software” on page 74
- “Set up a melt curve experiment in the software” on page 81
- “Set up a custom experiment in the software” on page 98

You can also use features of the desktop software to more easily set up some or all of an experiment. For example, you can set up an experiment using one of the following strategies.

- Use desktop software libraries to set up samples, targets or SNP assays, run methods, and analysis settings (see “Sample, target, and SNP assay libraries” on page 111).
- Import experiment parameters from external files or templates (see “Experiment setup using libraries, external files, and templates” on page 27).

Create or open a template

Create a new template or open an existing template in the  **Home** screen.

If security is enabled in the software, your account must have the permission of **Create New Template** to create a new template.

- In the  **New Experiment** pane, perform one of the following tasks to create a new template.

To	Action
Create a template without preexisting settings	Click Create New Experiment .
Create a template from a system template	- Select Create New Experiment ▶ Template. - Navigate to and select the desired file, then click Open. System template files are installed with the software in: <code><drive>:\Program Files (x86)\Applied BioSystems\QuantStudio Design & Analysis Software\templates</code>, where <i><drive></i> is the drive on which the software is installed.

- In the  **Open Existing Experiment** pane, perform the following tasks to open an existing template.
 - Click **Open**.
 - Navigate to and select the desired file, then click **Open**.

Enter template properties

If security is enabled in the software, your account must have the permission of **Edit Template Properties**.

If security is enabled in the software, your account must have the permission of **Edit Save As Destination** to select a different location to save the file.

- Click the **Properties** tab to open and edit the experiment properties.
- (Optional) In the **Name** field, modify the file name.

The **Name** field determines two file names:

- The initial EDT file name.

Note: After the initial EDT file save, modifying the **Name** field does not update the EDT file name.

To change the EDT file name *after* the initial save, click  **Save ▶ Save As**. In the **Save As** dialog box, enter the new EDT file name in the **File Name** field. Click **Browse** to edit the location to save the file, then click **Save**.

- The default file name for the EDS file created during an instrument run.
- (Optional) Click the **Barcode** field, then scan or enter a plate barcode.

4. (Optional) Enter information in the **User name** field and the **Comments** field, if applicable.
5. Select an **Instrument type**, **Block type**, **Experiment type**, **Chemistry** (reagents), and **Run Mode** (Fast or Standard cycling) from the dropdown lists.

Note: The experiment type defines the available options for the template setup. For more information on the parameters defined in each experiment type, see “Experiment definitions” on page 21 .

6. (Optional) Click **Manage chemistry details** (see “Enter reagent information”).
7. (Optional) Click  **Save** or click  **Save** ▶ **Save as**.
8. To save as, in the **Save As** dialog box, edit the file name in the **File Name** field, click **Browse** to edit the location to save the file, then click **Save**.

Enter reagent information

1. In the **Properties** tab, click **Manage chemistry details**.
2. Click  **Add**.
3. Enter the reagent type, name, part number, lot number, and expiration date.
4. (Optional) Add a custom attribute for a reagent.
 - a. Click  in the table header to add a column for a custom attribute.
 - b. Click the **Custom Attribute** column header, then enter a new attribute.
 - c. Select a cell in the **Custom Attribute** column, then enter its information.
 - d. (Optional) Click  in the header to delete a custom attribute from the table.
5. (Optional) Click  to delete a reagent from the table.
6. Click **Close**.

Scan a barcode using the optional barcode scanner

The instrument is compatible with an optional Handheld Barcode Scanner (Cat. No. [4488442](#), purchased separately). The barcode scanner reads Code 128 (alphanumeric), which supports 128 ASCII character barcodes.

1. Click the **Barcode** field.
2. Hold the scanner 20–30 cm away from a plate or container label and aim at the center of the barcode, then press the trigger.
3. Slowly move the scanning beam across the barcode until the scanner emits a high-pitched tone.

When the scanner scans a barcode, it automatically transmits the following information:

- Transmits the alphanumeric equivalent of the barcode to the barcode field.
- Transmits other reagent information (Lot #, Part #, Expiration Date, etc.)

For more information about the hand-held barcode scanner, see the user documentation provided with the barcode scanner.


Confirm or edit the run method and optical filter selection

In the **Method** tab, perform the follow tasks if needed.

If security is enabled in the software, your account must have the permission of **Edit Run Method** to edit the run method.

- (Optional) Adjust the reaction volume.
- (Optional) Edit the default run method (thermal protocol).
 - The default run method is optimized for TaqMan™ assays and a broad range of other reagents.
 - To edit the default run method, see “Adjust method parameters” on page 103.
- (Optional) Edit the default optical filter selection (see “Select optical filters” on page 105).
 - The default optical filter selection is for factory-calibrated (system) dyes.
 - For more information about instrument supported dyes and their calibration and optical filter selection, see the instrument user guide.

Method elements

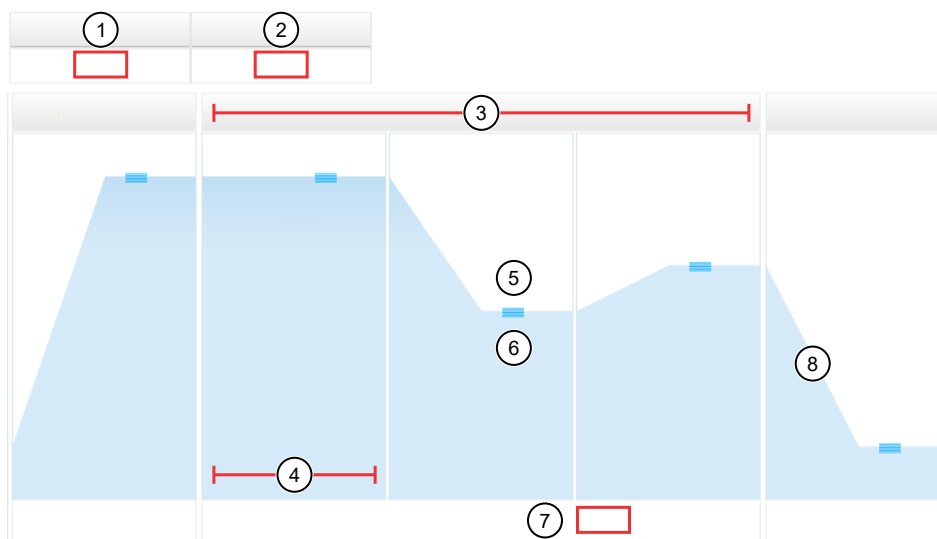

- ① Reaction volume
- ② Heated cover temperature
- ③ Stage of thermal protocol
- ④ Step within a stage

- ⑤ Step temperature
- ⑥ Step hold time
- ⑦ Number of cycles
- ⑧ Temperature ramp rate

Assign plate and well attributes using the Quick Setup subtab

If security is enabled in the software, your account must have the following permissions to perform the following tasks:

- The permission of **Assign Targets / Assays** to assign a target or an assay to a well
- The permission of **Define New Target / Assay** to add a new target or a new assay
- The permission of **Assign Sample** to assign a sample to a well
- The permission of **Define New Sample** to add a new sample
- The permission of **Edit Sample** to edit a sample
- The permission of **Delete Sample** to delete a sample
- The permission of **Edit Passive Reference** to edit the passive reference dye

Note: This section provides general procedures to set up the plate.

- For detailed procedures to set up a plate, see Appendix D, “Detailed procedures to set up plate / well details and libraries”.
 - For specific instructions for each experiment type, see the corresponding chapter in this guide.
-

1. In the **Plate** tab, select plate wells in the  **Plate Layout** or the  **Well Table**. For more information, see “Select plate wells” on page 21.
2. Click **Quick Setup**.
3. Assign the well attributes for the selected wells.
 - Enter the sample and target or SNP assay names in the text field.
 - Select a defined sample and target or SNP assay from the dropdown lists.
For more information about defining or importing samples and targets or SNP assays, see “Define and assign well attributes (Advanced Setup subtab)” on page 107.
 - Click **Advanced Setup**, then change the default selections for the reporter and quencher dyes and for tasks where applicable. For more information, see “Assign a task to wells” on page 108.
4. (Optional) Enter comments for the selected wells.
5. In the **Plate Attributes** pane, select a **Passive Reference** from the dropdown list.

Select plate wells

- Select plate wells in the  **Plate Layout**.

To	Action
Select a single well	Click a well in the plate
Select multiple wells	Click-drag in the plate
Select contiguous wells	Shift-click wells in the plate
Select non-contiguous wells	Ctrl-click wells in the plate
Select a column of wells	Click a column header
Select all wells	Click the top-left corner of the plate grid
Select a block of wells	Click a well to define a corner, then Shift-click another well on the opposite corner
Deselect a single well	Ctrl-click the selected well

- Select plate wells in the  **Well Table**.

To	Action
Select a single well	Click a row in the table
Select multiple wells	Click-drag in the table
Select contiguous wells	Shift-click rows in the table
Select non-contiguous wells	Ctrl-click rows in the table
Deselect a single well	Ctrl-click the selected row

Experiment definitions

The parameters that you can define vary by experiment type.

Experiment Type	Targets	SNP Assays	Samples	Biological Replicate Groups	Passive Reference	Reference and Endogenous Controls
Standard Curve	✓		✓	✓	✓	
Relative Standard Curve	✓		✓	✓	✓	✓
Comparative C _t	✓		✓	✓	✓	✓
Melt Curve	✓		✓		✓	
Genotyping		✓	✓		✓	
Presence / Absence	✓		✓		✓	
Custom	✓		✓	✓	✓	

Save a template file

Save a template file as an unlocked template

Note: You cannot save a locked template as an unlocked template.

If security is enabled in the software, your account must have the permission of **Edit Save As Destination** to select a different location to save the file.

- Save the template with the same EDT file name.
 - In any tab, click  **Save**.
 - In the menu bar, click **File** ▶ **Save**.
- Save the template with a new EDT file name.

Note: You cannot save a locked template with a new file name.

- In any tab, click  **Save** ▶ **Save as**.
- In the menu bar, click **File** ▶ **Save as**.

In the **Save As** dialog box, enter the new EDT file name in the **File Name** field. Click **Browse** to edit the location to save the file, then click **Save**.

Save a template file as a locked template

IMPORTANT! Always save a backup *unlocked* version of a template before saving it as a locked template.

Record the password for the locked template because lost passwords cannot be recovered.

If security is enabled in the software, your account must have the permission of **Edit Save As Destination** to select a different location to save the file.

1. In the menu bar, select **File** ▶ **Save As Locked Template** .
2. Enter and confirm a password, then click **OK**.

Note: The password is required to open the template with unlimited editing options. Without the password, a locked template can still be opened but with limited editing options.

3. In the **Save As** dialog box, edit the file name in the **File Name** field.
4. Click **Browse** to edit the location to save the file, then click **Save**.

Prepare reactions

See instrument user guide for information about compatible reagents and required materials for PCR reactions.

Follow the instructions provided by the manufacturer to prepare reactions.

Follow the other guidelines described in this section.

Good laboratory practices for PCR and RT-PCR


- Wear clean gloves and a clean lab coat.
 - Do not wear the same gloves and lab coat that you have previously used when handling amplified products or preparing samples.
- Change gloves if you suspect that they are contaminated.
- Maintain separate areas and dedicated equipment and supplies for:
 - Sample preparation and reaction setup.
 - Amplification and analysis of products.
- Do not bring amplified products into the reaction setup area.
- Open and close all sample tubes carefully. Avoid splashing or spraying samples.
- Keep reactions and components capped as much as possible.
- Use a positive-displacement pipettor or aerosol-resistant barrier pipette tips.
- Clean lab benches and equipment periodically with 10% bleach solution or DNA decontamination solution.

Guidelines for handling samples and reagents

- Use calibrated pipettors and aerosol-resistant tips.
- Prepare the reaction mixes according to the recommendations that are provided by the manufacturer of the master mixes and assay mixes.
- Include excess volume in calculations to account for the loss that occurs during reagent transfers.
- Use TE buffer or water to dilute samples and standards.
- Use care when diluting samples and standards. Mistakes or inaccuracies in making the dilutions affect data accuracy.
- Keep the dilutions and assay mix frozen and protected from light until use. Excessive exposure to light can affect the fluorescent probes or dyes.
- Perform the following tasks before each use:
 - Mix the master mix thoroughly by swirling the bottle.
 - Resuspend the assay mix by vortexing, then centrifuge the tube briefly.
 - Thaw any frozen samples by placing them on ice. When thawed, resuspend the samples by vortexing, then centrifuge the tubes briefly.

Guidelines for setting up the reactions in the plates or tubes

- Use good laboratory practices for PCR and RT-PCR. For more information, see “Good laboratory practices for PCR and RT-PCR” on page 23.
- Ensure that the arrangement of the PCR reactions matches the plate layout displayed in the software.
- Confirm that the liquid in each well is at the bottom of the well and free of bubbles. If it is not, centrifuge the plate again.

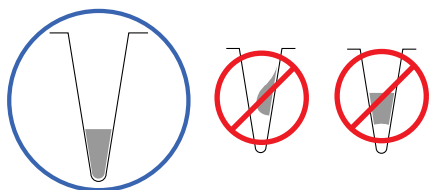

- Ensure that plates or tubes are properly sealed.
- Keep the reaction plate or tubes at 4°C and protected from light until you are ready to load the plate into the instrument.
- Keep the bottom of the plate clean. Fluids and other contaminants on the bottom of the plate can contaminate the sample block and cause an abnormally high background signal.
- If necessary, use a permanent marker or pen to mark a tube and the side of a plate. Do not use fluorescent markers.

Start and monitor a run

IMPORTANT! Before loading a reaction-filled plate into the instrument, review the detailed procedures (see “Load and unload the plate in the instrument” on page 94).

If security is enabled in the software, your account must have the permission of **Start Run**.

1. Go to an instrument connected to the computer that is running the desktop software.
2. Load the plate into the instrument.



CAUTION! (*QuantStudio™ 3 Real-Time PCR System and QuantStudio™ 5 Real-Time PCR System only*) The instrument should be used by trained operators who have been warned of the moving parts hazard.

3. Open a template (EDT file) in the desktop software.
4. (*Optional*) In the **Export** tab, select **Auto Export** to export run results automatically after the run ends.
5. In the **Method** tab and the **Plate** tab, review the template parameters and setup.
6. In the **Run** tab, select the instrument to use from the **START RUN** dropdown list.

7. Accept or edit the default name for the EDS file, then click **Save**.

Note: EDS files contain the run data and results. The system creates the EDS default file name for the EDS file from the **Experiment Name** in the **Properties** tab.

8. During the instrument run, monitor the run:

- In the **Run** tab of the desktop software.
- In the instrument touchscreen.

9. Unload the plate from the instrument.



CAUTION! PHYSICAL INJURY HAZARD. During instrument operation, the plate temperature can reach 100°C. Allow it to cool to room temperature before handling.

Note: (*QuantStudio™ 3 Real-Time PCR System and QuantStudio™ 5 Real-Time PCR System only*) If the instrument does not eject the plate, contact Support.

Note: If the connection between the instrument and the desktop software is interrupted during the run, the instrument still completes the run. However, the run data (EDS file) must be transferred from the instrument to the desktop software using a USB drive or a network drive.

View real-time run information in the desktop software

You can view real-time run information in the **Run** tab of the desktop software.

- In the **Run** tab, view the time remaining for the run and the run status.
- In the **Run** tab, in the **Amplification Plot** subtab, view real-time data and plots.
 - Click  **View** to select the data that are displayed in each well.
 - Select wells in the plate layout to highlight respective curves in the plot.
 - Select curves in the plot to highlight respective wells in the plate layout.

Note: For melt curve experiments, you can only monitor the melt curve plot.

View real-time run information on the instrument touchscreen

You can view real-time run information in the instrument home screen.

- In the instrument home screen, view the block temperature, the time remaining for the run, and the run status.
- Touch  or swipe left once to view real-time run method information.
- Touch  or swipe left twice to view real-time data and plots.

View real-time data and plots on the instrument touchscreen

1. In the instrument home screen, during an instrument run, touch ➤ or swipe left twice.
2. Touch **Well details**.
3. Touch **Samples**, **Targets**, or **Tasks** to select a graphical representation of each selection.
4. Touch **Close** to return to the home screen.

Adjust the display of real-time plots on the instrument touchscreen

1. In the instrument home screen, during an instrument run, touch ➤ or swipe left twice to view real-time data and plots.
2. Touch **Zoom**.
3. Touch 🔍 or 🔍 to zoom in or out.
4. Touch the arrows to pan left, right, up, or down on the graph.
5. Touch **Close** to return to the default view.

View the Post-Run Summary in the desktop software

You can view a summary of the run after the run ends.

In the **Run** tab, click the **Post-Run Summary** tab to view a summary of the run, including the following information:

- **Experiment Name**
- **User Name**
- **Errors Encountered**
- **Instrument Serial Number** and **Instrument Name**
- **Start Time**, **Stop Time**, and **Run Duration**

Experiment setup using libraries, external files, and templates

The desktop software offers the following features so that you can more easily set up some or all of an experiment using libraries, external files, and templates.

- Use desktop software libraries to set up samples, targets or SNP assays, run methods, and analysis settings (see “Libraries overview” on page 111).
- Import some or all of an experiment setup from external files or templates (see the following table).

Table 1 Experiment setup from external files or templates

Option	Action	Setup information
Import sample information (define samples).	Import a sample definition file (see “Assign samples using a sample definition file” on page 99).	– Sample name – <i>(Optional)</i> Custom sample properties
Import samples, targets, and well assignments.	Import a plate setup file (see “Assign samples and targets using plate setup files” on page 100).	Plate setup information: – Well number – Sample name – Sample color – Target name – Dyes – <i>(Optional)</i> Other well information
Set up the plate layout in a spreadsheet without saving to a special format. <i>or</i> Use a subset of the columns in a plate layout spreadsheet.	Copy-paste from an XLS file (see “Assign targets, samples, and biological replicate groups from an XLS file” on page 101).	Plate setup information: – Well number – Sample name – Biological Group – Target name – Task – Dyes – Quantity – Comments – <i>(Optional)</i> Other well information
Use a complete template setup from an existing EDT or EDS file.	Create a new template from an existing template or run results file (see “Create new EDT files using existing EDT and EDS files” on page 102).	– Plate setup information, as above – Reagent information – Thermal protocol – Analysis settings

3

General procedures to review results

■ About the quantification cycle (C_q)	28
■ Overview of the Results tab	29
■ Workflow: General procedures to review the run results	31
■ Guidelines for viewing and analyzing results	31
■ Assess results in the Amplification Plot	32
■ Assess results in the Well Table view	37
■ Review the dye signal profile using the Multicomponent Plot	38
■ Review the signal profile using the Raw Data Plot	39
■ Review the flags in the QC Summary	40
■ View calibration results in the desktop software	41
■ Export experiments or results	42

This section includes information about reviewing results and configuring analysis settings for all experiment types. For information about a specific experiment type, see the corresponding chapter in this guide.


For step-by-step instructions for results review procedures, see the desktop software Help.

About the quantification cycle (C_q)

Term	Name	Description
C_q	Quantification cycle	C_q is the general form for gene expression metrics. C_q values (both C_t and C_{rt}) are used as the primary input values for sample quantification experiments: absolute quantification (AQ) and relative quantification (RQ). C_t and C_{rt} are the algorithm-specific calculations of C_q.
C_t	Threshold cycle	The PCR cycle number at which the fluorescence signal meets the threshold in the amplification plot. C_t is the gene expression metric result when using the Baseline Threshold Algorithm.
C_{rt}	Relative threshold cycle	The PCR cycle number for the threshold calculated from the modeled amplification efficiency profile. C_{rt} is the gene expression metric result when using the Relative Threshold Algorithm.

Overview of the Results tab

Review and analyze run data in the **Results** tab. In the **Results** tab, two additional tools display at the right of the workflow bar.

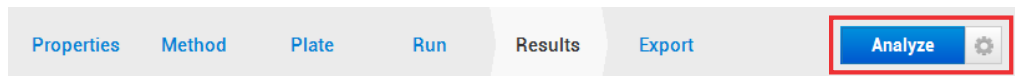

- Click **Analyze** after omitting wells or changing the analysis settings.
- Click to access analysis settings.

Note: The analysis settings and plots that are available vary by experiment type.

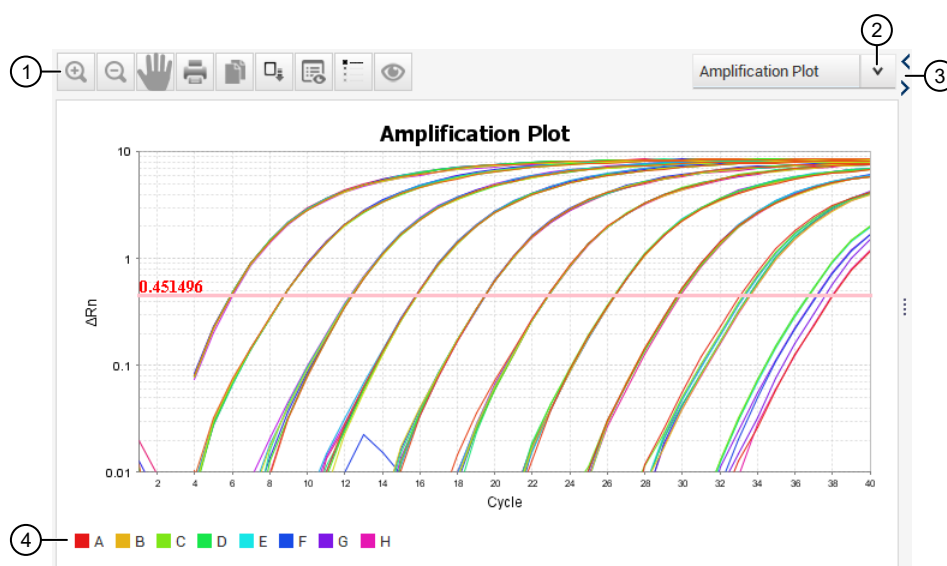


Figure 1 Plot pane

- | | |
|--|---|
| ① Plot toolbar | ③ Expand/contract the plot pane display |
| ② Plot selection list (varies by experiment) | ④ Plot legend |

The Plot toolbar includes the following options:

- | | |
|----------------------------|-----------------------------|
| • Zoom in and out | • Configure plot properties |
| • Print or copy plot image | • Show/hide plot legend |
| • Save plot as image file | • Configure plot settings |

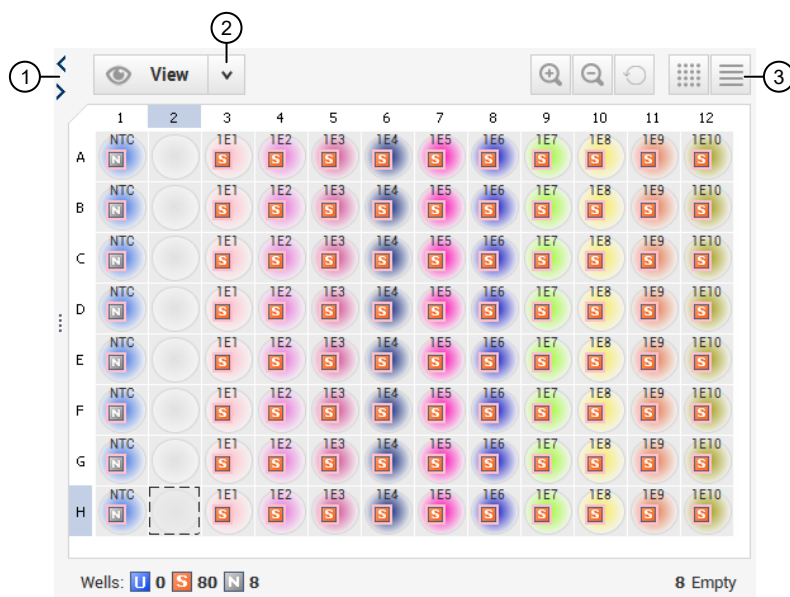


Figure 2 Plate Layout

- ① Expand/contract the **Plate Layout** display (layout is expanded in this figure)
- ② **View**: Select well properties to display
- ③ **Plate Layout** toolbar

The **Plate Layout** toolbar includes the following options:

- Zoom in and out
- Fit plate to window
- Display Plate Layout
- Display Well Table

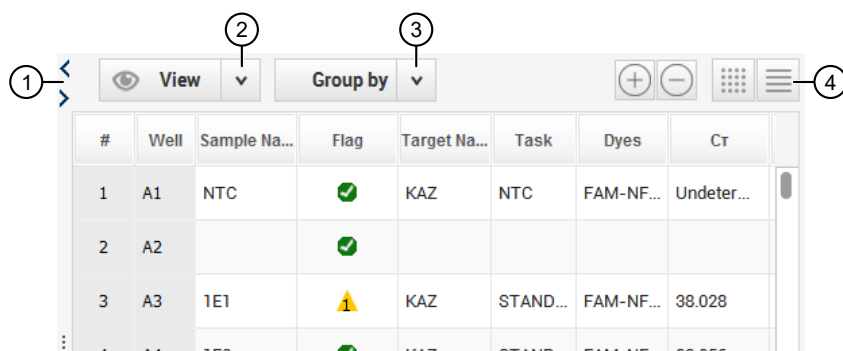


Figure 3 Well Table


- ① Expand/contract the **Well Table** display (table is expanded in this figure)
- ② **View**: Select well properties to display
- ③ **Group by**: Select a parameter by which to group well rows
- ④ **Well Table** toolbar

The **Well Table** toolbar includes the following options:

- Expand grouped rows
- Collapse grouped rows
- Display Plate Layout
- Display Well Table

Workflow: General procedures to review the run results

When a run is complete, the desktop software automatically analyzes the run data using the analysis settings that are specified during template development. The software then displays the run results in the **Results** tab.

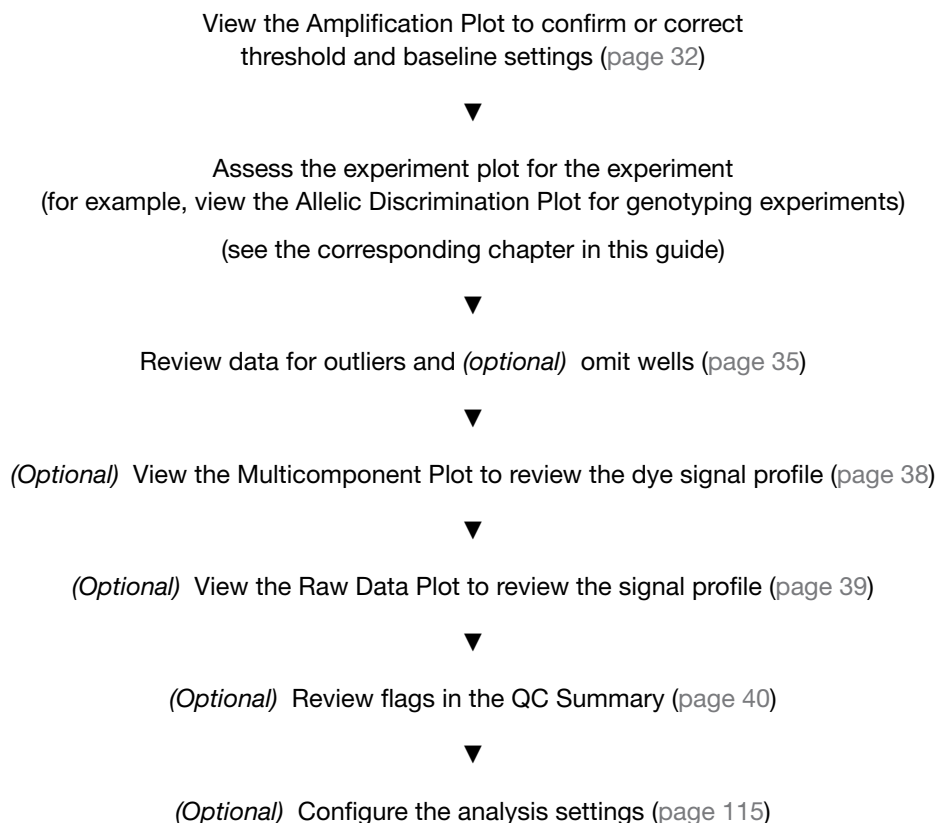

IMPORTANT! If you omit wells or configure the analysis settings, click **Analyze** to reanalyze the data.

Guidelines for viewing and analyzing results

- For information about adjusting the views in the **Results** tab, see the desktop software Help.
- To reanalyze the data, select all the wells in the  **Plate Layout**, then click **Analyze**.
- To enable auto-analysis of data after a run, select **Tools** ▶ **Preferences** ▶ **Experiment**, then select **Auto Analysis**.

Assess results in the Amplification Plot

Amplification Plot overview

The Amplification Plot displays sample amplification as a function of cycle number or well. You can use the amplification plot to perform the following tasks:

- Confirm or correct baseline and threshold values.
- Identify outliers.
- Identify and examine abnormal amplification. Abnormal amplification can exhibit one of the following characteristics:
 - Increased fluorescence in negative control wells.
 - Absence of detectable fluorescence at an expected cycle.

Note: If you notice abnormal amplification or a complete absence of fluorescence, see the instrument user guide for troubleshooting information.

Three plots are available. Some plots can be viewed as a linear or $\log_{10}$ graph.

Table 2 Amplification Plot types

Plot type	Description	Use to
ΔR_n vs Cycle	ΔR_n is the magnitude of normalized fluorescence signal, relative to the baseline fluorescence, generated by the reporter at each cycle during the PCR amplification.	• Identify and examine irregular amplification. • View threshold values for the run.
R_n vs Cycle	R_n is the fluorescence signal from the reporter dye normalized to the fluorescence signal from the passive reference.	• Identify and examine irregular amplification. • View baseline values for the run.
C_t vs Well	C_t is the PCR cycle number at which the fluorescence meets the threshold in the amplification plot.	• Locate outlying amplification (outliers).

Assess the overall shape of the Amplification Plot curves

You can assess the overall shape of the Amplification Plot curves in the **Results** tab.

If no data are displayed in the **Results** tab, click **Analyze**.

1. In the **Results** tab, select **Amplification Plot** from the dropdown list.
2. Click  to configure the plot, then make the following selections:
 - **Plot Type:** ΔR_n vs Cycle
 - **Graph Type:** Log
 - **Plot Color:** Target, Sample, or Well

The **Amplification Plot** is displayed for the selected wells in the  **Plate Layout**.

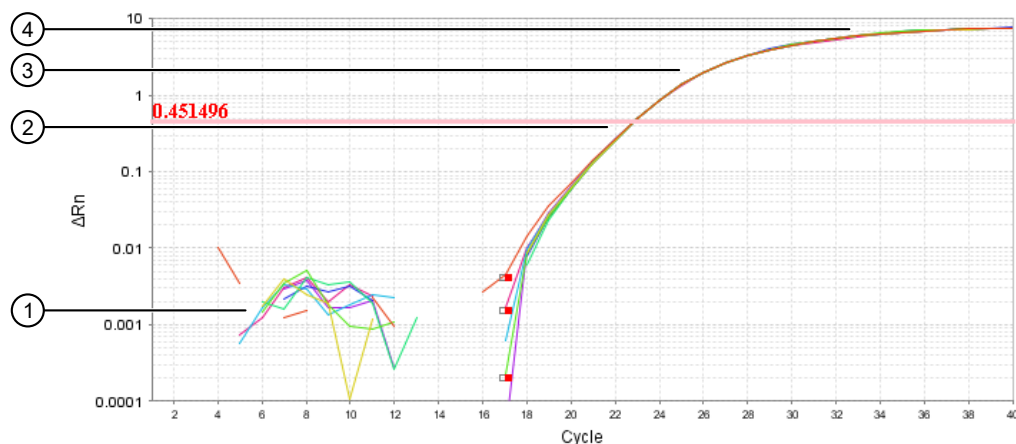


Figure 4 Typical Amplification Plot

A typical amplification curve has four distinct sections:

- | | |
|---------------------------------|-----------------|
| ① Baseline | ③ Linear phase |
| ② Exponential (geometric) phase | ④ Plateau phase |

Confirm or correct threshold settings




1. In the **Results** tab, select **Amplification Plot** from the dropdown list.
2. Click  to configure the plot, then make the following selections:
 - **Plot Type:** ΔR_n vs Cycle
 - **Graph Type:** Log
 - **Plot Color:** Target, Sample, or Well

The **Amplification Plot** is displayed for the selected wells in the  **Plate Layout**.

3. (Optional) Adjust the threshold.
 - Click-drag the threshold bar into the exponential phase of the curve.
 - Configure the C_t analysis settings (see “ C_t settings overview” on page 116).

Table 3 Examples of threshold settings

Set the threshold in the exponential phase of the amplification curve. A threshold set above or below the optimum can increase the standard deviation of the replicate groups.

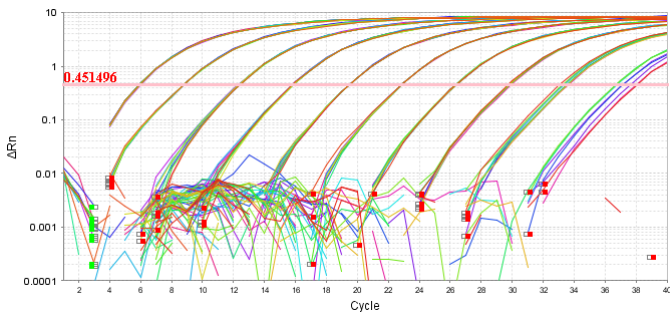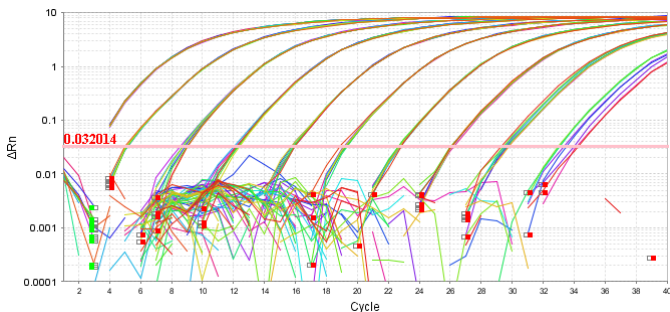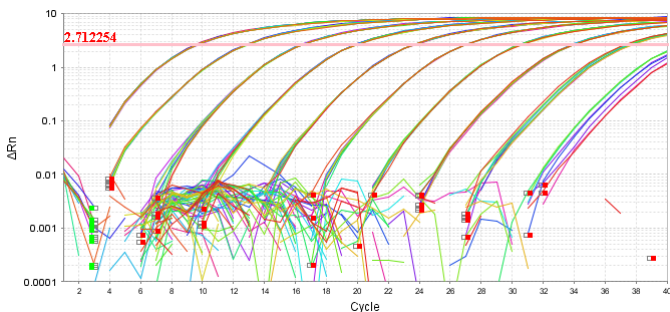Threshold setting evaluation	Example
Threshold set correctly.	
Threshold set too low.	
Threshold set too high.	

Confirm or correct baseline settings

1. In the **Results** tab, select **Amplification Plot** from the dropdown list.
2. Click  to configure the plot, then make the following selections:
 - **Plot Type:** Rn vs Cycle
 - **Graph Type:** Linear
 - **Plot Color:** Well
 - Select **Show: Baseline Start / Baseline End**

Note: The start and end cycles are used to calculate the baseline.

The **Amplification Plot** is displayed for the selected wells in the  **Plate Layout**.

The start (■) and end (■) cycles display for each well.

3. (Optional) Adjust the start and end cycle values for the baseline (see “Ct settings overview” on page 116).

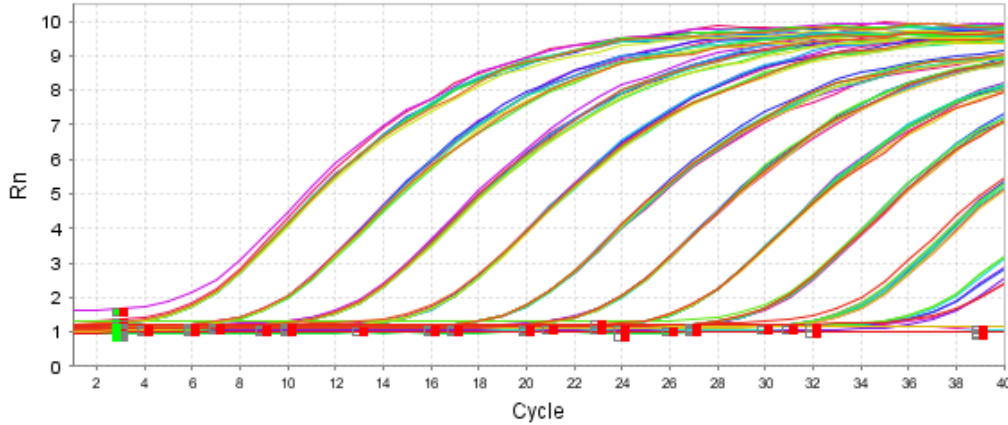


Figure 5 Example of correct baseline

Set the end cycle (■) a few cycles before the cycle number where significant fluorescence signal is detected.

Omit outliers from analysis

Outlier wells have C_q (C_t or C_{rt}) values that differ significantly from the average for the associated replicate wells. To ensure C_q (C_t or C_{rt}) precision, consider omitting the outliers from analysis.

If security is enabled in the software, your account must have the permission of **Omit Wells for Analysis**.

1. In the **Results** tab, select **Amplification Plot** from the dropdown list.
2. Click , then make the following selections to configure the plot:
 - **Plot Type: C_t vs Well**
 - **Graph Type: Linear**
 - **Plot Color: Well**

The C_t values are displayed for the selected wells in the  **Plate Layout**.

3. Click  to examine the **Well Table** for outliers.
 - a. Select **Group by ▶ Replicate**.
 - b. Identify outliers in each replicate group.
Outlier wells typically have one or more QC flags.

4. Omit outliers in either the  **Well Table** or  **Plate Layout** view.
 - In the  **Well Table**, select **Omit** in outlier rows of the table.
 - In the  **Plate Layout**, right-click a well, then select **Omit**.
5. Click **Analyze** to reanalyze the run data with any outliers removed.

Optimize display of negative controls in the Amplification Plot

1. In the **Results** tab, select **Amplification Plot** from the dropdown list.
2. Click  to configure the plot, then make the following selections:
 - **Plot Type: ΔR_n vs Cycle**
 - **Graph Type: Linear**
 - **Plot Color: Target**
 - Deselect **Show: Threshold**
 - Deselect **Show: Baseline Start / Baseline End**
3. In either the  **Plate Layout** or  **Well Table**, select the negative control wells (wells that should not have amplification for a particular target).
4. Click  (configure plot properties), then select the **Y Axis** tab.
 - a. Deselect **Auto-adjust range**.
 - b. Enter **Minimum value** of -1 .
 - c. Enter **Maximum value** of 2 .
 - d. Click **Save**.

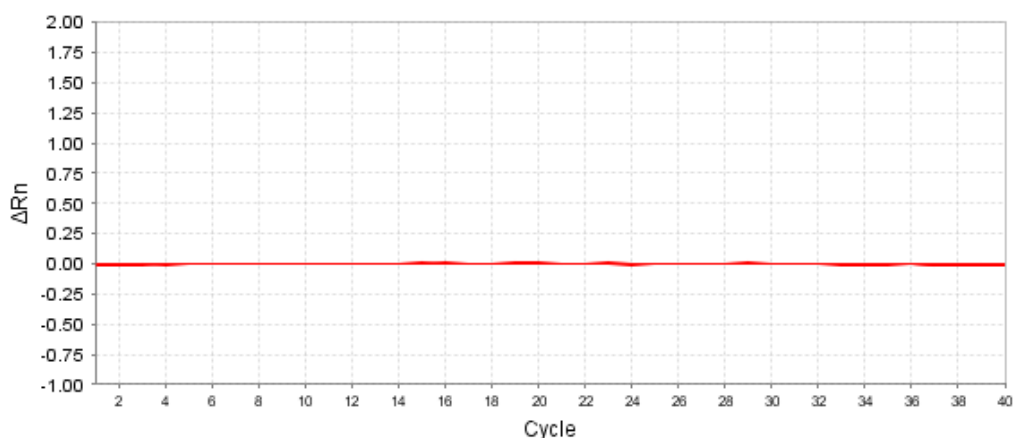


Figure 6 Example Amplification Plot of negative controls

The linear plot displays the Amplification Plot for negative controls as smooth lines. The expanded Y axis displays low levels of amplification.

Assess results in the Well Table view

Well Table overview

The  **Well Table** displays data for each well in the reaction plate. The data that are displayed depend on the specific experiment type and can include the following information:

- Sample name, target name, task, and dyes
- Values that are specific to particular stage of the method
For example: C_t or C_{rt} , normalized fluorescence (R_n), or melt temperature (T_m)
- Values that are specific to a particular experiment type
For example: genotype calls, presence/absence calls, or quantities
- Omitted wells
- QC flags
- Comments

Group or sort the Well Table

Some of the possible options for grouping or sorting the  **Well Table** are described in the following table. Available grouping categories depend on the specific experiment type and analysis settings.

If the data are grouped, this order is retained in the report (see “Print a report of run data and results” on page 44).

Note: You can select multiple columns when sorting, but you can only make one selection for grouping rows.

Group category	Description	Notes
Replicate ^[1,2,3]	Grouped by replicate	• Examine the C_t or quantity values for each replicate group to assess the precision of C_t values.
Flag	Grouped as flagged and unflagged wells	• A flag indicates that the software found a potential error in the flagged well. • For more information about QC flags, see the desktop software Help.
C_t ^[1,2,4]	Grouped by C_t value	• C_t value < 8—There may be too much template in the reaction. • C_t value > 35—There may be a low amount of target in the reaction; for C_t values > 35, expect a higher standard deviation.
RQ ^[2]	Grouped by RQ value	• RQ value < 1—There is less relative target in the test sample as compared to the calibrator sample. • RQ value > 1—There is more relative target in the test sample as compared to the calibrator sample.
Call ^[4,5]	Grouped by genotype call or presence/absence call	—

(continued)

Group category	Description	Notes
Tm1 ^[3]	Grouped by T _m value	Tm1 refers to the dominant peak. This grouping category is only applicable for the following experiment types: • Melt curve experiments • Any experiment with a melt curve data collection step (e.g., absolute standard curve)

^[1] For standard curve experiments.

^[2] For relative standard curve and comparative C_t ($\Delta\Delta C_t$) experiments.

^[3] For melt curve experiments.

^[4] For genotyping experiments.

^[5] For presence/absence experiments.

Review the dye signal profile using the Multicomponent Plot

Multicomponent Plot overview

The Multicomponent Plot displays the complete spectral contribution of each dye over the duration of the PCR run.

Use the Multicomponent Plot to obtain the following information.

- Confirm that the signal from the passive reference dye remains unchanged throughout the run.
- Review reporter dye signal for spikes, dips, and/or sudden changes.
- Confirm that no amplification occurs in the negative control wells.

View and assess the Multicomponent Plot

You can view and assess the **Multicomponent Plot** in the **Results** tab.

If no data are displayed in the **Results** tab, click **Analyze**.

1. In the **Results** tab, select **Multicomponent Plot** from the dropdown list.
2. Click  to configure the plot, then make the following selections:
 - **Plot Color: Dye**

The **Multicomponent Plot** is displayed for the selected wells in the  **Plate Layout**.

3. In the **Plate Layout**, select wells one at a time, then examine the **Multicomponent Plot** for the following plot characteristics.

Plot characteristic	Description
Passive reference dye	The passive reference dye fluorescence signal should remain relatively constant throughout the PCR process.
Reporter dye	The reporter dye fluorescence signal should display a flat region corresponding to the baseline, followed by a rapid rise in fluorescence as the amplification proceeds.
Irregularities in the signal	Spikes, dips, and/or sudden changes in the fluorescence signal may have an impact on the data.
Negative control wells	The negative control wells should show no significant increase in fluorescence signal.

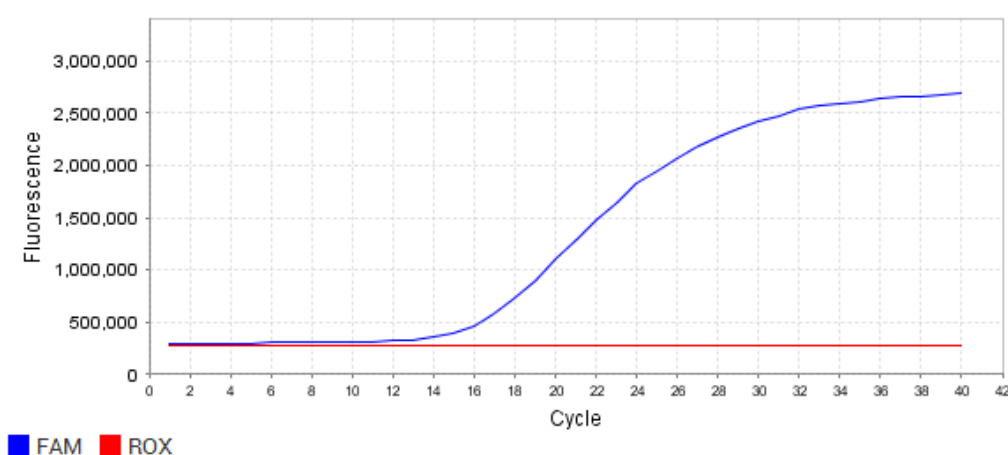


Figure 7 Example Multicomponent Plot (single well)

Review the signal profile using the Raw Data Plot

Raw Data Plot overview

The Raw Data Plot displays the raw fluorescence signal (not normalized) for each optical filter during each cycle of the real-time PCR.

View the Raw Data Plot to confirm a stable increase in signal (no abrupt changes or dips) from the appropriate filter.

View and assess the Raw Data Plot

You can view and assess the **Raw Data Plot** in the **Results** tab.

If no data are displayed in the **Results** tab, click **Analyze**.

1. In the **Results**, select **Raw Data Plot** from the dropdown list.
The **Raw Data Plot** is displayed for the selected wells in the **Plate Layout**.
2. Click  to display the **Show Cycle** scale.
3. Click-drag the **Show Cycle** pointer from cycle 1 to cycle 40, and confirm that each filter displays the characteristic signal increase.

For more information on each filter set, see the instrument user guide (see Appendix F, “Documentation and support”).



Figure 8 Example Raw Data Plot

Review the flags in the QC Summary

The **QC Summary** in the **Results** tab displays a list of the QC flags, including the flag frequency and location.

If no data are displayed in the **Results** tab, click **Analyze**.

1. In the **Results** tab, select **QC Summary** from the dropdown list.
2. Review the **Flag Details** table and the summary.
The **Wells** column of the **Flag Details** table identifies wells that triggered a flag.
3. (Optional) In the **Flag Details** table, click each flag to display a brief description of the flag.
For more information about the flag, see the desktop software Help.

View calibration results in the desktop software

Calibration results obtained on the instrument can be transferred via USB to the desktop software.

1. In the  **Home** screen, click  **Open**.
2. Navigate to and select the desired calibration EDS file.

If you are viewing calibration data files,  **Calibration QC Status** is displayed to the right of the  **Plate Layout** and  **Well Table**.

Note: For more information about using the  **Plate Layout** and  **Well Table**, see Chapter 3, “General procedures to review results” .

Review ROI/Uniformity calibration results

1. In the **ROI** tab, select a **Filter Set** from the dropdown list to see the corresponding results.
2. In the **Uniformity** tab, use the  **Plate Layout**,  **Well Table**, and  **Calibration QC Status** to review the run results.

Review Background calibration results

1. Select the plate wells in the  **Plate Layout** or the  **Well Table** to view the corresponding curves. For more information, see “Select plate wells” on page 21.
2. Review data in the  **Well Table**.
 - a. Review the results for each well in tabular format.
 - b. Sort the wells according to well or normalized fluorescence with each filter.
 - c. Select wells to review data in the analysis plot.
3. Click  **Calibration QC Status** to review the quality of the calibration data.

Review Dye calibration results

1. Select a Dye row in the **Calibration** table to view the corresponding analysis data plot.
2. Select the plate wells in the  **Plate Layout** or the  **Well Table** to view the corresponding curves. For more information, see “Select plate wells” on page 21.
3. Review data in the  **Well Table**.
 - a. Review the results for each well in tabular format.
 - b. Sort the wells according to well or normalized fluorescence with each filter.
 - c. Select wells to review data in the analysis plot.
4. Click  **Calibration QC Status** to review the quality of the calibration data.

Override the calibration data

Each EDS file contains the calibration data from the instrument on which it was run. You can use calibration data from another instrument for analysis of your run data.

- Calibration data must be from the same block type (96-well 0.2-mL block, 96-well 0.1-mL block, or 384-well block).
- Calibration data for a QuantStudio™ 1 Real-Time PCR Instrument run must be from a QuantStudio™ 1 Real-Time PCR Instrument run.
- Calibration data from a QuantStudio™ 5 Real-Time PCR Instrument run can override calibration data for a QuantStudio™ 3 Real-Time PCR Instrument run.
- Calibration data from a QuantStudio™ 3 Real-Time PCR Instrument run can *not* override calibration data for a QuantStudio™ 5 Real-Time PCR Instrument run.

1. Open the EDS file to recalibrate.
2. Select **Analysis ▶ Override Calibration ▶ Use Calibration From Another File**.
3. Navigate to, then select the EDS file containing the alternative calibration data.
4. Click **Open**.
5. (Optional) To revert to the original calibration data, select **Analysis ▶ Override Calibration ▶ Revert To Original Calibration**.

Export experiments or results

For step-by-step instructions for exporting experiments or results, see the desktop software Help.

Options for exporting the plots

To	Action
Save a plot as an image file	In the Results tab, click  (Save to image file). The  (Save to image file) button is located above the plot. Use the dropdown list above the plot to select a plot type.
Print a plot	Click  (Print). The  (Print) button is located above the plot. Use the dropdown list above the plot to select a plot type.
Copy a plot to the clipboard	Click  (Copy). The  (Copy) button is located above the plot. Use the dropdown list above the plot to select a plot type.

Export run data results

You can customize the exported results. See “Customize exported results” on page 43.

If security is enabled in the software, your account must have the permission of **Edit Export Destination** to edit the location in the **Export Location** field.

1. Select the **Export** tab.
2. (Optional) Edit the file name in the **File Name** field.
The default file name is set up in the **Name** field of the **Properties** tab.

3. Select a file type from the **File Type** dropdown list.

- **QuantStudio**
- **RDML**

For more information, see “Export configurations” on page 45.

4. If **QuantStudio** was selected as the file type, select a file format from the dropdown list.
For more information, see “Export configurations” on page 45.

5. Select a location to save the file in the **Export Location** field.

6. Select the **Open exported files when complete** checkbox if you would like to open the file immediately.

7. In the **Content** pane, select the checkboxes that are associated with the information to include in the export.

Selecting the content is not available if **RDML** was selected as the file type.

8. In the **Options** pane, select the radio button to include all of the data in a single file or to export separate files.

- **Unify the above content into one file**
- **Split the above content items into individual files**

For a XLS file, the content is separated onto different tabs if one file is selected.

For a TXT file, the content is included sequentially in the file if one file is selected.

The file names are appended with the content type if individual files are selected.

9. Click **Export**.

Customize exported results

Ensure that the content is selected. The options that are selected determine what is customized.

1. In the **Export tab**, click **Customize**.
The **Customize** dialog box is displayed. Each tab in the dialog box corresponds to content item that was selected.
2. In the **Customize** dialog box, select the tab that corresponds to the content to customize.

3. In the **Select Content** pane, select items to include in the exported file.
4. In the **Results** tab, select the following checkboxes to skip the corresponding data from being exported.
 - **Skip Empty Wells** checkbox
 - **Skip Omitted Wells** checkbox
5. Click **Close** to return to the **Export** tab.

Print the plate layout

1. From any tab, click **File ▶ Print**.
2. In the **Print Plate Layout** dialog box, select the radio button that corresponds with the color.
 - **Use sample color** radio button
 - **Use task color** radio button
3. Click **Print** to open the system dialog box for the printer.
4. Click **Preview** to open the **Print Preview** dialog box.
A preview of the plate layout is displayed.
5. In the **Print Preview** dialog box, click one of the following options.
 - Click **Save As PDF** to save the file.
 - Click **Print Report** to open the system dialog box for the printer.

Send run results to PowerPoint™

1. Open an EDS file.
2. In the menu bar, click **File ▶ Send to PowerPoint**.
3. In the **Create Slides** dialog box, select the content to include in the slides, then click **Apply**.

Note: Analysis plots are not available for slides until the run status is complete and the run data are analyzed.

Print a report of run data and results

If security is enabled in the software, your account must have the permission of **Edit Report Destination** to select a different location to save the file.

1. Open an EDS file.
2. In the menu bar, select **File ▶ Print Report**.
3. Select the run data and results to include in the report.
The available data and results vary according to the experiment type.

If the data were grouped by a specific parameter, that order is retained in the report. For information about grouping, see Figure 3 on page 30 and “Group or sort the Well Table” on page 37.

4. Click **Print Report** to open the system dialog box for the printer.
5. (Optional) Click **Print Preview**.
A preview of the report is displayed in the **Print Preview** dialog box.
6. In the **Print Preview** dialog box, click one of the following options.
 - Click **Save As PDF** to save the file.
 - Click **Print Report** to open the system dialog box for the printer.
7. In the **Save As** dialog box, edit the file name in the **File Name** field.
8. Click **Browse** to edit the location to save the file, then click **Save**.

Export configurations

Data type	Description	File format
Plate setup files, for future experiments	Plate setup information For example, the well number, sample name and color, target name, dyes, and other reaction plate contents.	• XLS • XLSX • TXT
Analyzed data, for further analysis	QuantStudio™ format	• XLS • XLSX • TXT
	RDML (Real-Time PCR Data Markup Language) format Used for standard curve, relative standard curve, and comparative C _t experiments.	• RDML

4

Set up, run, and review standard curve experiments

■ Standard curve experiments	46
■ Set up a standard curve experiment in the software	47
■ Set up and run the PCR reactions	48
■ Review results	49

Standard curve experiments

Overview

Use standard curve experiments to determine absolute target quantity in samples.

In a standard curve experiment, the software performs the following tasks.

1. The software measures amplification of the target in a standard dilution series and in test samples.
2. The software generates a standard curve using data from the standard dilution series.
3. The software uses the standard curve to interpolate the absolute quantity of target in the test samples.

Reaction types

Multiple targets can be assayed in a standard curve experiment, but each target requires its own standard curve.

Table 4 Reaction types for standard curve experiments

Reaction type (task)	Sample description
Standard	A sample that contains known or known relative quantities of the target • For known quantities—Quantify the target in the standard sample using an independent method. • For known relative quantities—Generate a relative dilution series of the target standards.

Table 4 Reaction types for standard curve experiments *(continued)*

Reaction type (task)	Sample description
Unknown	Test sample
No-template control (NTC/ Negative Control)	Water or buffer No amplification of the target should occur in NTC wells.

- The precision of quantification experiments improves as the number of replicate reactions increases. Set up the number of replicates appropriate for the experiment.
- For accurate and precise efficiency measurements, set up the standard dilution series with at least five dilution points over a broad range of standard quantities, 4 to 6 logs (10^4 - to 10^6 -fold). A concentrated template, such as a plasmid or PCR product, is best for this purpose.
A narrow range of standard quantities may be appropriate if the amount of standard is limited, the target is in low abundance, or the target is known to fall within a given range.

Compatible PCR options

Table 5 PCR options for standard curve experiments

Single- or multiplex PCR	PCR or RT-PCR ^[1]	Detection chemistry
Singleplex	PCR	TaqMan™
Multiplex	1-step RT-PCR 2-step RT-PCR	SYBR™ Green

^[1] RT-PCR: reverse transcription-PCR

Set up a standard curve experiment in the software

1. In the  **Home** screen, create or open a template.
 - In the  **New Experiment** pane, click **Create New Experiment** to create a new template.
 - In the  **Open Existing Experiment** pane, click **Open** to select and open an existing template.
2. In the **Properties** tab, enter the template information.
3. In the **Method** tab, adjust the reaction volume.
For most experiments, the default run method is appropriate.
4. In the **Plate** tab (**Quick Setup**), assign plate attributes.
 - a. In the **Plate Attributes** pane, select the **Passive Reference** from the dropdown list.
5. In the **Plate** tab (**Quick Setup**), define and assign well attributes.
 - a. Select wells in the  **Plate Layout** or the  **Well Table**.

b. Assign samples and targets to selected wells.

- Enter new sample and target names in the text fields.
- Select previously defined samples and targets from the dropdown lists.

Note: New sample or target names entered in the **Quick Setup** subtab are automatically populated with default values for **Reporter** (FAM) and **Quencher** (NFQ-MGB) and assigned a task (**U Unknown**). Edit these values in the **Advanced Setup** subtab.

6. (Optional) In the **Plate** tab, set up standard dilutions (see “Define and set up standard dilutions” on page 109).

7. (Optional) In the **Plate** tab (**Advanced Setup**), assign tasks.

a. Select wells in the  **Plate Layout** or the  **Well Table**.

b. In the **Targets** table, select the checkbox of a target, then select a task from the dropdown list.

Reaction type	Task
Unknown (test sample)	
No-template control	

8. (Optional) In the **Plate** tab (**Advanced Setup**), define and assign biological replicate groups (see “Define and assign biological replicate groups” on page 110).

In the **Plate** tab (**Advanced Setup**), ensure the **Samples** table contains the following sample information:


- One sample name for each technical replicate group of an unknown sample
- (Optional) A standard sample for each target

Set up and run the PCR reactions

1. Assemble the PCR reactions. Follow the manufacturer's instructions for the reagents and follow the plate layout set up in the software (see “Prepare reactions” on page 23).
2. In the desktop software, open the appropriate template (EDT file).
3. Load the reaction plate into the instrument and start the run (see “Start and monitor a run” on page 24).

Review results

Workflow: Review standard curve experiments

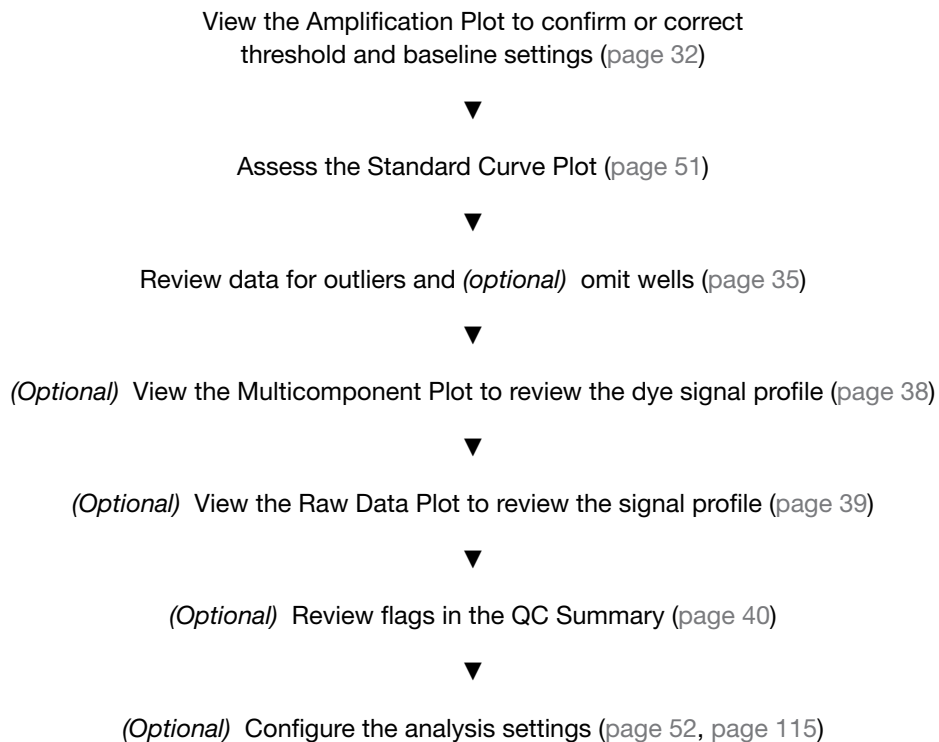

IMPORTANT! If you omit wells or configure the analysis settings, click **Analyze** to reanalyze the data.

Standard Curve Plot overview

The Standard Curve Plot displays the standard curve for samples designated as standards. The software calculates the quantity of an unknown target from the standard curve.

Table 6 Results or metrics to review in the Standard Curve Plot

Results or metrics	Description	Criteria for evaluation
Slope and amplification efficiency	The amplification efficiency is calculated using the slope of the regression line in the standard curve.	A slope close to -3.3 indicates optimal, 100% PCR amplification efficiency. Factors that affect amplification efficiency: • Improper design of the primer and probe • Range of standard quantities—For accurate and precise efficiency measurements, use a broad range of standard quantities, 5 to 6 logs (10^5- to 10^6-fold). • Number of standard replicates—For accurate efficiency measurements, include replicates to decrease the effects of pipetting inaccuracies. • PCR inhibitors—PCR inhibitors and contamination in the reaction can reduce amplification efficiency. • Other possible factors: – Component and properties of the reaction mix, such as salt content, DMSO, pH, etc. – Inaccurate sample or reagent pipetting – Improper analysis settings – Incorrect plate setup
R^2 value (correlation coefficient)	The R^2 value is a measure of the closeness of fit between the regression line and the individual C_q data points of the standard reactions.	• A value of 1.00 indicates a perfect fit between the regression line and the data points. • An R^2 value > 0.99 is desirable.
Error	The standard error of the slope of the regression line in the standard curve. The error can be used to calculate a confidence interval (CI) for the slope and therefore the amplification efficiency.	Acceptable value is determined by the experimental criteria.

Table 6 Results or metrics to review in the Standard Curve Plot (continued)

Results or metrics	Description	Criteria for evaluation
C _t values	The threshold cycle (C _t) is the PCR cycle number at which the fluorescence level meets the threshold.	A C_t value > 8 and < 35 is desirable. - C_t value < 8—There may be too much template in the reaction. - C_t value > 35—There may be a low amount of target in the reaction; for C_t values > 35, expect a higher standard deviation.

View and assess the Standard Curve Plot

You can view and assess the Standard Curve Plot in the **Results** tab.

If no data are displayed in the **Results** tab, click **Analyze**.

1. Select **Standard Curve** from the dropdown list.
2. Click  to configure the plot, then make the following selections:
 - **Target:** Select the target of interest
 - **Plot Color:** **Sample**, **Target**, or **Task**
 - Select all wells in the  **Plate Layout**

The Standard Curve Plot is displayed. The slope, R² value, amplification efficiency, and error are displayed below the plot.

3. Confirm that the slope, R² value, amplification efficiency, and error meet the experimental criteria.
4. Visually check that all unknown sample C_q (C_t or C_{rt}) values fall within the standard curve range.
5. In the  **Well Table**, use the **Group By** dropdown list to confirm that the C_q (C_t or C_{rt}) values of all replicate samples meet the experimental criteria.

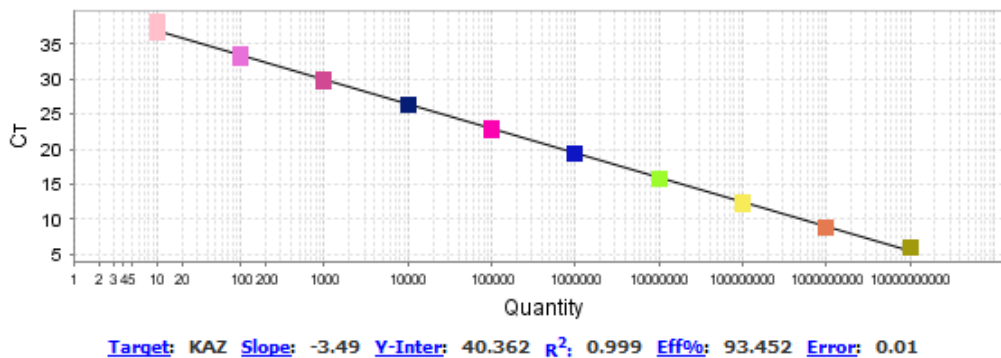


Figure 9 Example Standard Curve Plot

If the results do not meet the experimental criteria, troubleshoot using one of the following strategies:

- Omit wells, then reanalyze. For more information, see “Omit outliers from analysis” on page 35.
- Repeat the experiment, adjusting the template setup and analysis settings to improve results.

To learn more about the Standard Curve Plot, see “Standard Curve Plot overview” on page 50.

Standard curve settings overview

You can use the standard curve from another experiment and apply it to the current experiment. The two experiments must be from the same instrument type, block type, and run method.

To import an external standard curve, select  **Analysis Settings** ▶ **Standard Curve Settings**, then follow the instructions on the screen.

For step-by-step instructions for adjusting the standard curve settings, see the desktop software Help.



Set up, run, and review relative standard curve experiments and comparative C_t experiments

■ Relative standard curve experiments	53
■ Comparative Ct experiments	55
■ Relative quantitation: relative standard curve vs. comparative Ct	56
■ Set up a relative standard curve experiment in the software	57
■ Set up a comparative Ct experiment in the software	58
■ Set up and run the PCR reactions	59
■ Review results	60

Relative standard curve experiments

Overview

Use relative standard curve experiments to determine relative target quantity in samples.

In a relative standard curve experiment, the software performs the following tasks.

1. The software measures amplification of the target of interest and of an endogenous control target in a standard dilution series, in a reference (calibrator) sample, and in test samples.
The endogenous control is a target that is expressed equally in all samples; examples of endogenous controls are β -actin, GAPDH, and 18S ribosomal RNA. The software can algorithmically incorporate multiple endogenous control targets in relative quantification calculations.
The reference sample is used as the basis for relative quantification results (or 1 $\times$ sample). For example, in a study of drug effects on gene expression, an untreated control is an appropriate reference sample.
2. The software generates standard curves for the target of interest and the endogenous control using data from the corresponding standard dilution series.
3. The software uses the standard curves to interpolate the quantities of the target of interest and the endogenous control in each sample. The target quantity in each sample is then normalized to the sample's endogenous control quantity.
4. To determine the relative quantity of the target in test samples, the software divides the normalized target quantity in the sample by the normalized target quantity in the reference sample.

For a comparison of this method to the comparative C_t ($\Delta\Delta C_t$) method, see “Relative quantitation: relative standard curve vs. comparative C_t ” on page 56.

Reaction types

Relative standard curve experiments include the following reaction types for the endogenous control target and each target of interest.

Table 7 Reaction types for relative standard curve experiments

Reaction type (task)	Sample description
Standard	A sample that contains known or known relative quantities of the target - For known quantities—Quantify the target in the standard sample using an independent method. - For known relative quantities—Generate a relative dilution series of the target standards.
Reference sample ^[1]	The sample that is used as the basis for relative quantification results
Unknown	Test or reference sample
No-template control (NTC/ Negative Control)	Water or buffer No amplification of the target should occur in NTC wells.

^[1] To identify a sample as a reference sample, review the relative quantification settings.

- The precision of quantification experiments improves as the number of replicate reactions increases. Set up the number of replicates appropriate for the experiment.
- For accurate and precise efficiency measurements, set up the standard dilution series with at least five dilution points over a broad range of standard quantities, 4 to 6 logs (10^4 - to 10^6 -fold). A concentrated template, such as a plasmid or PCR product, is best for this purpose.
A narrow range of standard quantities may be appropriate if the amount of standard is limited, the target is in low abundance, or the target is known to fall within a given range.

Compatible PCR options

Table 8 PCR options for relative standard curve experiments

Single- or multiplex PCR	PCR or RT-PCR ^[1]	Detection chemistry
Singleplex	PCR	TaqMan™
Multiplex	1-step RT-PCR	SYBR™ Green
	2-step RT-PCR	

^[1] RT-PCR: reverse transcription PCR.

Comparative C_t experiments

Overview

Use comparative C_t ($\Delta\Delta C_t$) experiments to determine relative target quantity in samples.

In a comparative C_t experiment, the software performs the following tasks.

1. The software measures amplification of the target of interest and of an endogenous control target in a reference (calibrator) sample and in test samples.

The endogenous control is a target that is expressed equally in all samples; examples of endogenous controls are β -actin, GAPDH, and 18S ribosomal RNA. The software can algorithmically incorporate multiple endogenous control targets in relative quantification calculations.

The reference sample is used as the basis for relative quantification results (or $1\times$ sample). For example, in a study of drug effects on gene expression, an untreated control is an appropriate reference sample.

2. The measurements for the target of interest are normalized to the endogenous control.
3. To determine the relative quantity of the target in test samples, the software compares the normalized ΔC_q (ΔC_t or ΔC_{rt}) for the sample to the normalized ΔC_q (ΔC_t or ΔC_{rt}) for the reference sample.

For a comparison of this method to the relative standard curve method, see “Relative quantitation: relative standard curve vs. comparative C_t ” on page 56.

Reaction types

Comparative C_t experiments include the following reaction types for the endogenous control target and each target of interest.

Table 9 Reaction types for comparative C_t experiments

Reaction type (task)	Sample description
Reference sample ^[1]	The sample that is used as the basis for relative quantification results
Unknown	Test or reference sample
No-template control (NTC/ Negative Control)	Water or buffer No amplification of the target should occur in NTC wells.

^[1] To identify a sample as a reference sample, review the relative quantification settings.

The precision of quantification experiments improves as the number of replicate reactions increases. Set up the number of replicates appropriate for the experiment.

Compatible PCR options

Table 10 PCR options for comparative C_t experiments

Single- or multiplex PCR	PCR or RT-PCR ^[1]	Detection chemistry
Singleplex	PCR	TaqMan™
Multiplex	1-step RT-PCR	SYBR™ Green
	2-step RT-PCR	

^[1] RT-PCR: reverse transcription PCR

Relative quantitation: relative standard curve vs. comparative C_t

Use either relative standard curve experiments or comparative C_t experiments to determine the relative quantity of a target of interest in a test sample relative to a reference sample. Relative quantitation experiments are commonly used for the following applications.

- Comparison of expression levels of a gene in different tissues.
- Comparison of expression levels of a gene in a treated sample vs. an untreated sample.
- Comparison of expression levels of a gene of interest in different genetic backgrounds.
- Analysis of the gene expression changes over time under specific treatment conditions.

Table 11 Comparison of relative standard curve experiments and comparative C_t experiments

Characteristic	Relative standard curve	Comparative C_t
Typical use	Best for assays that have suboptimal PCR efficiency.	Best for high-throughput measurements of relative gene expression of many genes in many samples.
Advantage	Requires the least amount of validation because the PCR efficiencies of the target and endogenous control do not need to be equivalent.	• Relative levels of target in samples can be determined without the use of a standard curve, if the PCR efficiencies of the target and endogenous control are relatively equivalent. • Reduced reagent usage. • More space available in the reaction plate.
Limitation	A standard curve must be constructed for each target, which requires more reagents and more space in the reaction plate.	• Suboptimal (low PCR efficiency) assays may produce inaccurate results. • Before you use the comparative C_t method, we recommend that you determine that the PCR efficiencies for the target assay and the endogenous control assay are approximately equal.

Set up a relative standard curve experiment in the software

1. In the  **Home** screen, create or open a template.
 - In the  **New Experiment** pane, click **Create New Experiment** to create a new template.
 - In the  **Open Existing Experiment** pane, click **Open** to select and open an existing template.
2. In the **Properties** tab, enter the template information.
3. In the **Method** tab, adjust the reaction volume.
For most experiments, the default run method is appropriate.
4. In the **Plate** tab (**Quick Setup**), assign plate attributes.
 - a. In the **Plate Attributes** pane, select a **Passive Reference**, **Reference Sample**, and **Endogenous Control** from the dropdown lists.
5. In the **Plate** tab (**Quick Setup**), define and assign well attributes.
 - a. Select wells in the  **Plate Layout** or the  **Well Table**.
 - b. Assign samples and targets to selected wells.
 - Enter new sample and target names in the text fields.
 - Select previously defined samples and targets from the dropdown lists.

Note: New sample or target names entered in the **Quick Setup** subtab are automatically populated with default values for **Reporter** (FAM) and **Quencher** (NFQ-MGB) and assigned a task ( **Unknown**). Edit these values in the **Advanced Setup** subtab.

6. In the **Plate** tab, set up standard dilutions (see “Define and set up standard dilutions” on page 109).
7. (Optional) In the **Plate** tab (**Advanced Setup**), assign tasks.
 - a. Select wells in the  **Plate Layout** or the  **Well Table**.
 - b. In the **Targets** table, select the checkbox of a target, then select a task from the dropdown list.

Reaction type	Task
Unknown (test sample)	
No-template control	

8. (Optional) In the **Plate** tab (**Advanced Setup**), define and assign biological replicate groups (see “Define and assign biological replicate groups” on page 110).

The targets of interest and the endogenous control target should each have wells that are assigned with standard dilutions, unknown, and no-template-control tasks, and corresponding samples.

In the **Plate** tab (**Advanced Setup**), ensure the **Samples** table contains the following samples.

- Unknown samples
- Reference sample
- (Optional) For a dilution of target standards, each dilution step for each endogenous control target and target of interest has its own sample name.

Set up a comparative C_t experiment in the software

1. In the  **Home** screen, create or open a template.
 - In the  **New Experiment** pane, click **Create New Experiment** to create a new template.
 - In the  **Open Existing Experiment** pane, click **Open** to select and open an existing template.
2. In the **Properties** tab, enter the template information.
3. In the **Method** tab, adjust the reaction volume.
For most experiments, the default run method is appropriate.
4. In the **Plate** tab (**Quick Setup**), assign plate attributes.
 - a. In the **Plate Attributes** pane, select a **Passive Reference**, **Reference Sample**, and **Endogenous Control** from the dropdown lists.
5. In the **Plate** tab (**Quick Setup**), define and assign well attributes.
 - a. Select wells in the  **Plate Layout** or the  **Well Table**.
 - b. Assign samples and targets to selected wells.
 - Enter new sample and target names in the text fields.
 - Select previously defined samples and targets from the dropdown lists.

Note: New sample or target names entered in the **Quick Setup** subtab are automatically populated with default values for **Reporter** (FAM) and **Quencher** (NFQ-MGB) and assigned a task ( **Unknown**). Edit these values in the **Advanced Setup** subtab.

6. (Optional) In the **Plate** tab (**Advanced Setup**), assign tasks.
 - a. Select wells in the  **Plate Layout** or the  **Well Table**.
 - b. In the **Targets** table, select the checkbox of a target, then select a task from the dropdown list.

Reaction type	Task
Unknown (test sample)	
No-template control	

7. (Optional) In the **Plate** tab (**Advanced Setup**), define and assign biological replicate groups (see “Define and assign biological replicate groups” on page 110).

The targets of interest and the endogenous control target should each have wells assigned with unknown and no-template-control tasks, and corresponding samples.

In the **Plate** tab (**Advanced Setup**), ensure the **Samples** table contains:

- Unknown samples
- Reference sample

Set up and run the PCR reactions

1. Assemble the PCR reactions. Follow the manufacturer's instructions for the reagents and follow the plate layout set up in the software (see “Prepare reactions” on page 23).
2. In the desktop software, open the appropriate template (EDT file).
3. Load the reaction plate into the instrument and start the run (see “Start and monitor a run” on page 24).

Review results

Workflow: Review relative standard curve and comparative C_t experiments

Relative standard curve experiments

Comparative C_t experiments

View the Amplification Plot to confirm or correct threshold and baseline settings (page 32)



Assess the Standard Curve Plot (page 51)



Review data for outliers and *(optional)* omit wells (page 35)



Assess the Gene Expression Plot (page 63)



(Optional) View the Endogenous Control Profile using the QC Plot (page 64)



(Optional) View the Multicomponent Plot to review the dye signal profile (page 38)



(Optional) View the Raw Data Plot to review the signal profile (page 39)



(Optional) Review flags in the QC Summary (page 40)



(Optional) Configure the analysis settings (page 65, page 115)

IMPORTANT! If you omit wells or configure the analysis settings, click **Analyze** to reanalyze the data.

View and assess the Standard Curve Plot

This section only applies to relative standard curve experiments.

Standard Curve Plot overview

The Standard Curve Plot displays the standard curve for samples designated as standards. The software calculates the quantity of an unknown target from the standard curve.

Table 12 Results or metrics to review in the Standard Curve Plot

Results or metrics	Description	Criteria for evaluation
Slope and amplification efficiency	The amplification efficiency is calculated using the slope of the regression line in the standard curve.	A slope close to -3.3 indicates optimal, 100% PCR amplification efficiency. Factors that affect amplification efficiency: • Improper design of the primer and probe • Range of standard quantities—For accurate and precise efficiency measurements, use a broad range of standard quantities, 5 to 6 logs (10^5- to 10^6-fold). • Number of standard replicates—For accurate efficiency measurements, include replicates to decrease the effects of pipetting inaccuracies. • PCR inhibitors—PCR inhibitors and contamination in the reaction can reduce amplification efficiency. • Other possible factors: – Component and properties of the reaction mix, such as salt content, DMSO, pH, etc. – Inaccurate sample or reagent pipetting – Improper analysis settings – Incorrect plate setup
R^2 value (correlation coefficient)	The R^2 value is a measure of the closeness of fit between the regression line and the individual C_q data points of the standard reactions.	• A value of 1.00 indicates a perfect fit between the regression line and the data points. • An R^2 value > 0.99 is desirable.
Error	The standard error of the slope of the regression line in the standard curve. The error can be used to calculate a confidence interval (CI) for the slope and therefore the amplification efficiency.	Acceptable value is determined by the experimental criteria.

Table 12 Results or metrics to review in the Standard Curve Plot (continued)

Results or metrics	Description	Criteria for evaluation
C_t values	The threshold cycle (C_t) is the PCR cycle number at which the fluorescence level meets the threshold.	A C_t value > 8 and < 35 is desirable. C_t value < 8—There may be too much template in the reaction. C_t value > 35—There may be a low amount of target in the reaction; for C_t values > 35, expect a higher standard deviation.

View and assess the Standard Curve Plot

You can view and assess the Standard Curve Plot in the **Results** tab.

If no data are displayed in the **Results** tab, click **Analyze**.

1. Select **Standard Curve** from the dropdown list.
2. Click  to configure the plot, then make the following selections:
 - **Target:** Select the target of interest
 - **Plot Color:** **Sample**, **Target**, or **Task**
 - Select all wells in the  **Plate Layout**

The Standard Curve Plot is displayed. The slope, R^2 value, amplification efficiency, and error are displayed below the plot.

3. Confirm that the slope, R^2 value, amplification efficiency, and error meet the experimental criteria.
4. Visually check that all unknown sample C_q (C_t or C_{rt}) values fall within the standard curve range.
5. In the  **Well Table**, use the **Group By** dropdown list to confirm that the C_q (C_t or C_{rt}) values of all replicate samples meet the experimental criteria.

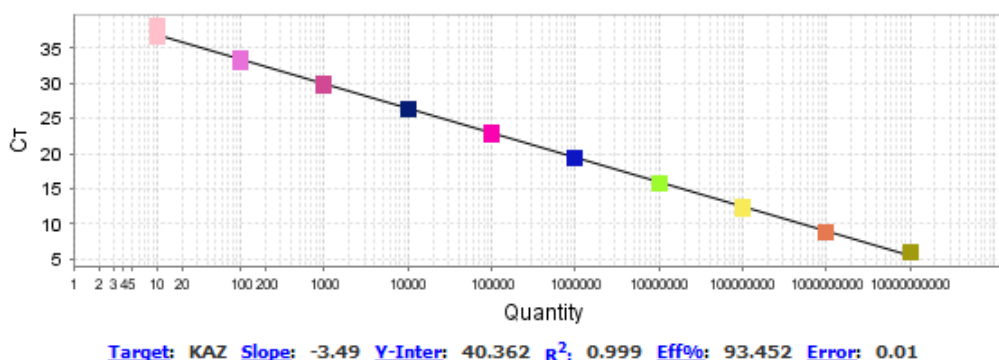


Figure 10 Example Standard Curve Plot

If the results do not meet the experimental criteria, troubleshoot using one of the following strategies:

- Omit wells, then reanalyze. For more information, see “Omit outliers from analysis” on page 35.
- Repeat the experiment, adjusting the template setup and analysis settings to improve results.

To learn more about the Standard Curve Plot, see “Standard Curve Plot overview” on page 50.

Standard curve settings overview

You can use the standard curve from another experiment and apply it to the current experiment. The two experiments must be from the same instrument type, block type, and run method.

To import an external standard curve, select **Analysis Settings** ▶ **Standard Curve Settings**, then follow the instructions on the screen.

For step-by-step instructions for adjusting the standard curve settings, see the desktop software Help.

Gene Expression Plot overview

The Gene Expression Plot displays the results of relative quantification calculations for relative standard curve and comparative C_t experiments.

Review the Gene Expression Plot to evaluate the fold change in expression level of the targets of interest in the test samples relative to the reference sample.

There are two plots available, depending on the experimental focus. Each plot can be viewed on a linear, $\log_{10}$, $\ln$, and $\log_2$ scale.

Table 13 Gene Expression plots

Plot type	Description
RQ vs. Target	Groups the relative quantification (RQ) values by target. Each sample is plotted for each target.
RQ vs. Sample	Groups the relative quantification (RQ) values by sample. Each target is plotted for each sample.

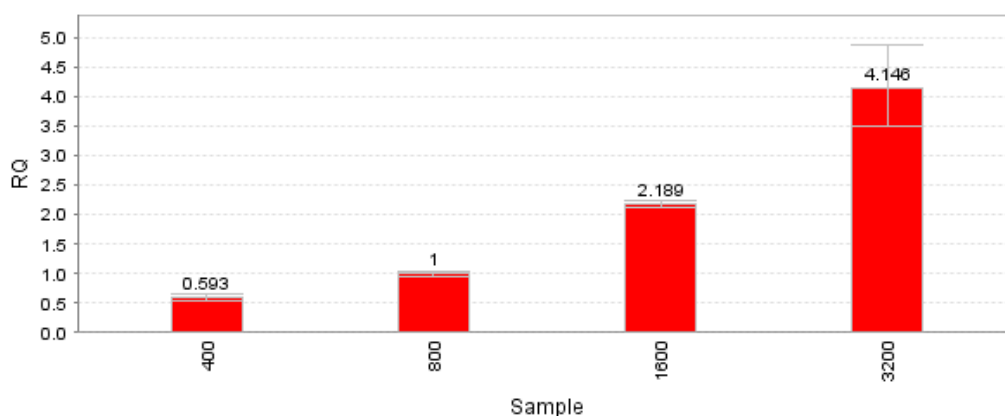


Figure 11 Example Gene Expression Plot

In this example, there is one target of interest, and the reference sample (calibrator) is sample 800.

QC Plot overview

The QC Plot is a visual display of the C_q (C_t or C_{rt}) levels of potential endogenous control targets across all samples (Endogenous Control Profile).

Use the QC Plot to help choose the best endogenous control for an experiment. Select the target with a quantity (indicated by C_q (C_t or C_{rt}) value) that does not change under experimental conditions.

All targets can be displayed in the QC Plot. You can view up to four potential endogenous controls at a time.

View and assess the QC Plot

You can view and assess the **QC Plot** in the **Results** tab.

If no data are displayed, click **Analyze**.

1. In the **Results** tab, select **QC Plot** from the dropdown list.
The **Endogenous Control Profile** is displayed.
2. In the right pane, select the targets to display, then select the color and shape from the dropdown lists.
3. (Optional) In the **View Replicate Results Table** tab, select the samples to omit from analysis.
4. (Optional) To change the endogenous controls used for analysis, select **Analysis Settings** ▶ **Relative Quantification Settings** (see “Relative quantification settings overview” on page 65).
5. Click **Analyze** to see the result of the adjustments.

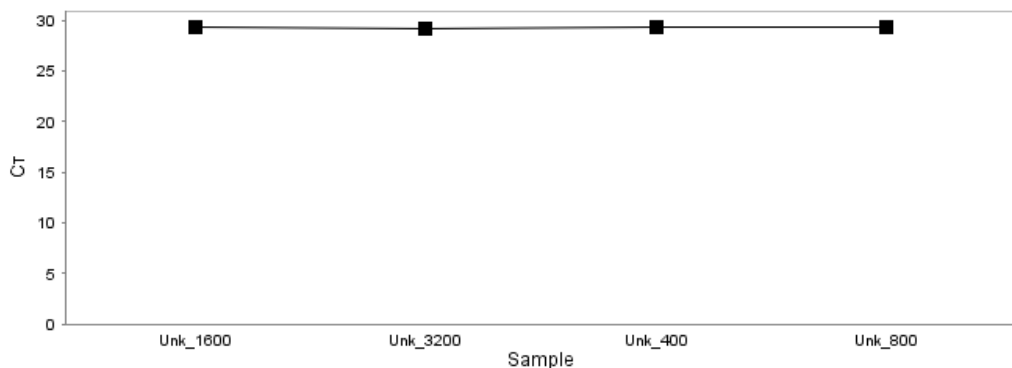


Figure 12 Example QC Plot

Relative quantification settings overview

In the **Results** tab, select  **Analysis Settings** ▶ **Relative Quantification Settings** to configure the following parameters:

Parameter	Description
Analysis Type	Select Multiplex or Singleplex analysis.
References	Set the reference sample, or set a biological replicate group as the reference sample.
Endogenous Controls	Change the endogenous control, or select multiple endogenous controls.
Efficiency (Comparative C_t experiments only)	Set the amplification efficiency for a target. The amplification efficiency for each target is calculated from the standard dilution series in relative standard curve experiments.
Outlier Rejection (Multiplex reactions only)	Outliers with ΔC_q (ΔC_t or ΔC_{rt}) values less than or equal to the entered value are rejected.
RQ Min/Max Calculations	Determines the algorithm used to calculate the relative quantification minimum and maximum values (error bars). • Confidence Level—Select to calculate the RQ minimum and maximum values based on the selected confidence level. • Standard Deviations—Select to calculate the RQ minimum and maximum values based on the selected number of standard deviations.

For step-by-step instructions for adjusting the relative quantification settings, see the desktop software Help.



Set up, run, and review genotyping experiments

■ Genotyping experiments	66
■ Set up a genotyping experiment in the software	67
■ Set up and run the PCR reactions	68
■ Review results	69

Genotyping experiments

Overview

Use genotyping experiments to detect single nucleotide polymorphism (SNP) variants of a target nucleic acid sequence.

Genotyping experiments use preformulated TaqMan™ SNP Genotyping Assays that include the following components:

- Two sequence-specific primers for amplification of sequences containing the SNP of interest
- Two allele-specific TaqMan™ probes for Allele 1 and Allele 2

In a genotyping experiment, the software performs the following tasks.

1. The software normalizes the fluorescence of the reporter dyes to the fluorescence of the passive reference dye in each well.
2. The software plots the normalized reporter dye signal of each sample well on an Allelic Discrimination Plot, which contrasts the reporter dye intensities of the allele-specific probes.
3. The software algorithmically clusters the sample data, and assigns a genotype call to the samples of each cluster according to its position on the plot.

Reaction types

Table 14 Reaction types for genotyping experiments

Reaction type (task)	Sample description
Unknown	Test sample
No-template control	Water or buffer No amplification of the target should occur in NTC wells.
Allele control (1/1)	Control sample that is homozygous for allele 1
Allele control (1/2)	Control sample that is heterozygous allele 1/allele 2
Allele control (2/2)	Control sample that is homozygous for allele 2

Allele controls are optional but recommended. Including allele controls helps to improve the clustering algorithm, particularly in situations where a limited number of samples are run.

In genotyping experiments, the software makes calls for individual wells. Running 3 or more replicates of each reaction can help identify outlier wells that may be present.

Compatible PCR options

Table 15 PCR options for genotyping experiments

Single- or multiplex PCR	PCR or RT-PCR ^[1]	Detection chemistry
Multiplex ^[2]	PCR	TaqMan™

^[1] RT-PCR: reverse transcription PCR

^[2] Each SNP genotyping assay is a multiplex assay with a probe for each allele. Multiple SNP assays can be performed in a single well.

Genotyping calls are based either on end-point data (data collected outside of any PCR cycling stage) or on real-time data (data collected during a PCR cycling stage). For detailed information, see the analysis settings section of this guide.

We recommend collecting real-time amplification data during the PCR stage, for troubleshooting purposes.

Set up a genotyping experiment in the software

1. In the  **Home** screen, create or open a template.
 - In the  **New Experiment** pane, click **Create New Experiment** to create a new template.
 - In the  **Open Existing Experiment** pane, click **Open** to select and open an existing template.
2. In the **Properties** tab, enter the template information.
3. In the **Method** tab, adjust the reaction volume.
For most experiments, the default run method is appropriate.

4. In the **Plate** tab (**Quick Setup**), assign plate attributes.
 - a. In the **Plate Attributes** pane, select the **Passive Reference** from the dropdown list.
5. In the **Plate** tab (**Quick Setup**), define and assign well attributes.
 - a. Select wells in the  **Plate Layout** or the  **Well Table**.
 - b. Assign samples and SNP assays to selected wells.
 - Enter new sample and SNP assay names in the text fields.
 - Select previously defined samples and SNP assays from the dropdown lists.

Note: New sample or SNP assay names entered in the **Quick Setup** subtab are automatically populated with the following default values:

Allele 1	Reporter	Quencher	Allele 2	Reporter	Quencher	Task
Allele 1	VIC	NFQ-MGB	Allele 2	FAM	NFQ-MGB	 Unknown

Edit these values in the **Advanced Setup** subtab.

6. (Optional) In the **Plate** tab (**Advanced Setup**), assign tasks.
 - a. Select wells in the  **Plate Layout** or the  **Well Table**.
 - b. In the **SNP Assays** table, select the checkbox of a SNP assay, then select a task from the dropdown list.

Reaction type	Task
Unknown (test sample)	
No-template control	
Allele control (1/1) ^[1]	
Allele control (1/2) ^[1]	
Allele control (2/2) ^[1]	

^[1] Optional but recommended

In the **Plate** tab (**Advanced Setup**), ensure the **Samples** table contains the following samples:

- Unknown samples
- (Optional) Allele control samples

Set up and run the PCR reactions

1. Assemble the PCR reactions. Follow the manufacturer's instructions for the reagents and follow the plate layout set up in the software (see “Prepare reactions” on page 23).
2. In the desktop software, open the appropriate template (EDT file).
3. Load the reaction plate into the instrument and start the run (see “Start and monitor a run” on page 24).

Review results

Workflow: Review genotyping experiments



IMPORTANT! If you omit wells or configure the analysis settings, click **Analyze** to reanalyze the data.

Allelic Discrimination Plot overview

The Allelic Discrimination Plot contrasts the R_n or the ΔR_n of the reporter dyes for the allele-specific probes of the SNP assay. It is an intermediary step in the software algorithm for genotyping calls.

Data points tend to cluster along the horizontal axis (Allele 1), vertical axis (Allele 2), or diagonal (Allele 1/Allele 2).

Table 16 Data clusters in the Allelic Discrimination Plot

A substantial increase in...	Clusters along...	Indicates...
Fluorescence of VIC™ dye-labeled probe only	Horizontal axis	Homozygosity for Allele 1
Fluorescence of FAM™ dye-labeled probe only	Vertical axis	Homozygosity for Allele 2
Fluorescence of both VIC™ and FAM™ dye-labeled probes	Diagonal	Heterozygosity for Allele 1 – Allele 2

Review the allelic discrimination plot to assess data clusters.

- Confirm that clustering of control samples is as expected.
- Visually assess clusters for the three possible genotypes.

Note: The desktop software clustering algorithm does not call genotypes if all the samples are one genotype (form one cluster).

View and assess the Allelic Discrimination Plot

You can view and assess the **Allelic Discrimination Plot** in the **Results** tab.

If no data are displayed in the **Results** tab, click **Analyze**.

1. In the **Results** tab, select **Allelic Discrimination Plot** from the dropdown list.
2. Click  to configure the plot, then make the following selections:
 - **SNP Assay:** select the assay of interest
 - **Plot Type:** **Cartesian** or **Polar**

The Allelic Discrimination Plot is displayed for the selected SNP assay.

Note: Initially, all points in the plot are cyan because all of the wells in the  **Plate Layout** are selected. Click anywhere in the plot or  **Plate Layout** to deselect all wells. The data points in the plot change to the call colors.

3. Confirm that control data clusters as expected.
 - a. In the  **Well Table** or  **Plate Layout**, select the wells containing a control to highlight the corresponding data points in the plot.
 - b. Check that the data points for each genotype control cluster along the expected axis of the plot.
4. Select the cluster at the bottom-left corner of the plot, then confirm that only the negative control wells are selected in the  **Plate Layout** or  **Well Table**.
 Samples can unexpectedly cluster with the negative controls for one of the following reasons.
 - Samples contain no DNA.
 - Samples contain PCR inhibitors.
 - Samples are homozygous for a sequence deletion.
5. Review the other clusters in the plot.
 - a. Click-drag a box around a cluster to select the associated wells.
 - b. Confirm that the expected wells are selected in the  **Plate Layout** or  **Well Table**.
6. Look for outliers outside the three genotype clusters.

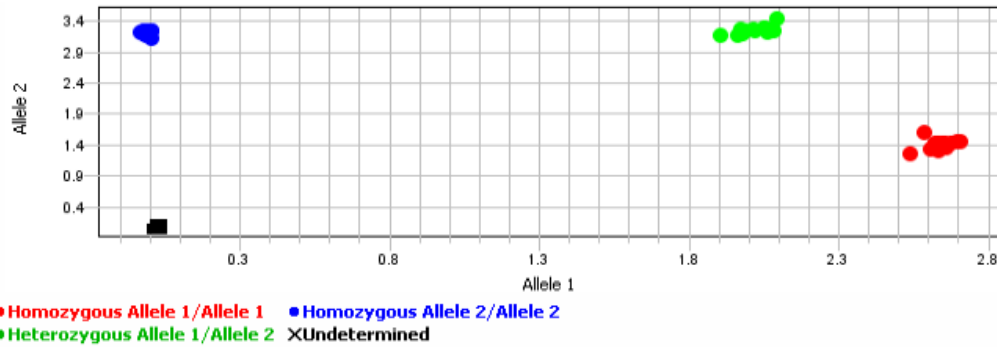


Figure 13 Example Allelic Discrimination Plot

To confirm results, retest outliers and samples with no amplification (cluster with negative controls).

Perform manual calls

You can perform manual calls in the **Results** tab.

1. In the **Results** tab, select **Allelic Discrimination Plot** from the dropdown list.
2. If the data are not analyzed, click **Analyze**.
3. (For multiple assays only) Click , then select a SNP assay from the dropdown list.
4. In the **Allelic Discrimination Plot**, use the lasso tool to select the samples to be manually called.
5. Click , then select the allele call from the **Apply Call** dropdown list.
6. Click **Analyze**.

IMPORTANT! To maintain manual calls after reanalysis, select  **Analysis Settings** ▶ **Call Settings**, then deselect **Default Settings** and select **Keep Manual Calls from Previous Analysis**.

Note: To remove manual calls, select  **Analysis Settings** ▶ **Call Settings**, deselect **Keep Manual Calls from Previous Analysis**, then reanalyze.

Call settings overview (genotyping)

In the **Results** tab, select  **Analysis Settings** ▶ **Call Settings** to edit the following settings:

- Data analysis settings
- Default call settings for SNP assays without custom call settings
- Custom call settings for individual SNP assays

Table 17 Options for data analysis settings (genotyping experiments)

Data analysis setting	Description
Analyze Data from Post-PCR Read Only	Only post-PCR read data is used to determine calls.
Analyze Data from Pre-PCR Read and Post-PCR Read ^[1]	The pre-PCR read is subtracted from the post-PCR read to determine calls.
Analyze Real-Time Rn Data ^[2,3]	The normalized reporter data (Rn) from the user-selected cycle of the cycling stage is used to determine calls.
Analyze Real-Time Rn - Median (Rna to Rnb) ^[2,4,3]	A <i>quick baseline-subtracted Rn</i> from the user-selected cycle of the cycling stage is used to determine calls. The quick baseline-subtracted Rn is the Rn minus the median value of the baseline region. The median subtraction provides improved data accuracy.
Analyze Real-Time dRn Data ^[2,3]	The regular ΔRn (dRn) from the user-selected cycle of the cycling stage is used to determine calls. The ΔRn is calculated by subtracting the best-fit line through the baseline region. This method is better if the baselines are not flat.

^[1] The run method must include a pre-read stage.

^[2] Data collection must be on during the PCR stage.

^[3] Analysis is not restricted to the last cycle; adjust the analysis cycle using the **Reveal Traces** feature while viewing the Allelic Discrimination Plot.

^[4] Rna to Rnb refers to all the cycles from the Start Cycle Number to the End Cycle Number of the baseline region.

Table 18 Options for call settings (genotyping experiments)

Call setting	Description
Autocaller Enabled	The autocaller algorithm is used to make genotype calls.
Keep Manual Calls from Previous Analysis	If autocaller is enabled, maintains manual calls after reanalysis.
Quality Value	The Quality Value is a proprietary estimation of the likelihood that a genotyping call is correct (associated with the correct cluster). If the Quality Value is less than the setting, the call is undetermined.

For step-by-step instructions for adjusting the call settings, see the desktop software Help.

7

Set up, run, and review presence/absence experiments

■ Presence/absence experiments	73
■ Set up the presence/absence experiment in the software	74
■ Set up and run the PCR reactions	75
■ Review results	76

Presence/absence experiments

Overview

Use presence/absence experiments to determine the presence or absence of a target nucleic acid sequence in a sample.

The software calls the target present or absent based on an algorithmically determined call threshold. (The call threshold is different from the C_t threshold; the C_t threshold is not used to make calls.)

Reaction types

Presence/absence reaction types depend on whether the experiment is set up with or without an internal positive control (IPC).

- **Presence/absence experiments with IPC (recommended)** are multiplex assays for the target of interest and the IPC target. The IPC is used to confirm that a negative result for the target of interest is not caused by a failed PCR.

Table 19 Reaction types for presence/absence experiments with IPC

Reaction type (task)	Sample description
Unknown	Test sample <i>and</i> IPC template
Negative control	Water or buffer <i>and</i> IPC template
No amplification control (NAC; blocked IPC) ^[1]	Water or buffer plus a blocking agent <i>and</i> IPC template; amplification prevented by blocking agent

^[1] Minimum of two replicates is required for this control.

- **Presence/absence experiments without IPC** are singleplex reactions.

Table 20 Reaction types for presence/absence experiments without IPC

Reaction type (task)	Sample description
Unknown	Test sample
Negative control	Water or buffer

The software makes calls for individual wells. Running three or more replicates of each reaction can help identify outlier wells that may be present.

Compatible PCR options

Table 21 PCR options for presence/absence experiments

Single- or multiplex PCR	PCR or RT-PCR ^[1]	Detection chemistry
Singleplex (without IPC)	PCR	TaqMan™
Multiplex (with IPC)	1-step RT-PCR 2-step RT-PCR	

^[1] RT-PCR: reverse transcription-PCR

Presence/absence calls are based on end-point data (data collected after the PCR stage).

- The data collected is the normalized intensity of the reporter dye, or Rn.
- If end-point experiments include pre-PCR data points, the software calculates the delta Rn (ΔRn) value according to the following formula:

$$\Delta Rn = Rn_{(post-PCR \text{ read})} - Rn_{(pre-PCR \text{ read})}$$
 , where Rn = normalized readings.

We recommend collecting real-time amplification data during the PCR stage, for troubleshooting purposes.

Set up the presence/absence experiment in the software

1. In the  **Home** screen, create or open a template.
 - In the  **New Experiment** pane, click **Create New Experiment** to create a new template.
 - In the  **Open Existing Experiment** pane, click **Open** to select and open an existing template.
2. In the **Properties** tab, enter the template information.
3. In the **Method** tab, adjust the reaction volume.
For most experiments, the default run method is appropriate.
4. In the **Plate** tab (**Quick Setup**), assign plate attributes.
 - a. In the **Plate Attributes** pane, select the **Passive Reference** from the dropdown list.

5. In the **Plate** tab (**Quick Setup**), define and assign well attributes.
 - a. Select wells in the  **Plate Layout** or the  **Well Table**.
 - b. Assign samples and targets to selected wells.
 - Enter new sample and target names in the text fields.
 - Select previously defined samples and targets from the dropdown lists.

Note: New sample or target names entered in the **Quick Setup** subtab are automatically populated with default values for **Reporter** (FAM) and **Quencher** (NFQ-MGB) and assigned a task ( **Unknown**). Edit these values in the **Advanced Setup** subtab.

6. (Optional) In the **Plate** tab (**Advanced Setup**), assign tasks.
 - a. Select wells in the  **Plate Layout** or the  **Well Table**.
 - b. In the **Targets** table, select the checkbox of a target, then select a task from the dropdown list.


Reaction type	Target	Task
Unknown (test sample)	Target of interest	
	IPC	
Negative control	Target of interest	
	IPC	
NAC (blocked IPC)	Target of interest	
	IPC	

Set up and run the PCR reactions

1. Assemble the PCR reactions. Follow the manufacturer's instructions for the reagents and follow the plate layout set up in the software (see “Prepare reactions” on page 23).
2. In the desktop software, open the appropriate template (EDT file).
3. Load the reaction plate into the instrument and start the run (see “Start and monitor a run” on page 24).

Review results

Workflow: Review presence/absence experiments

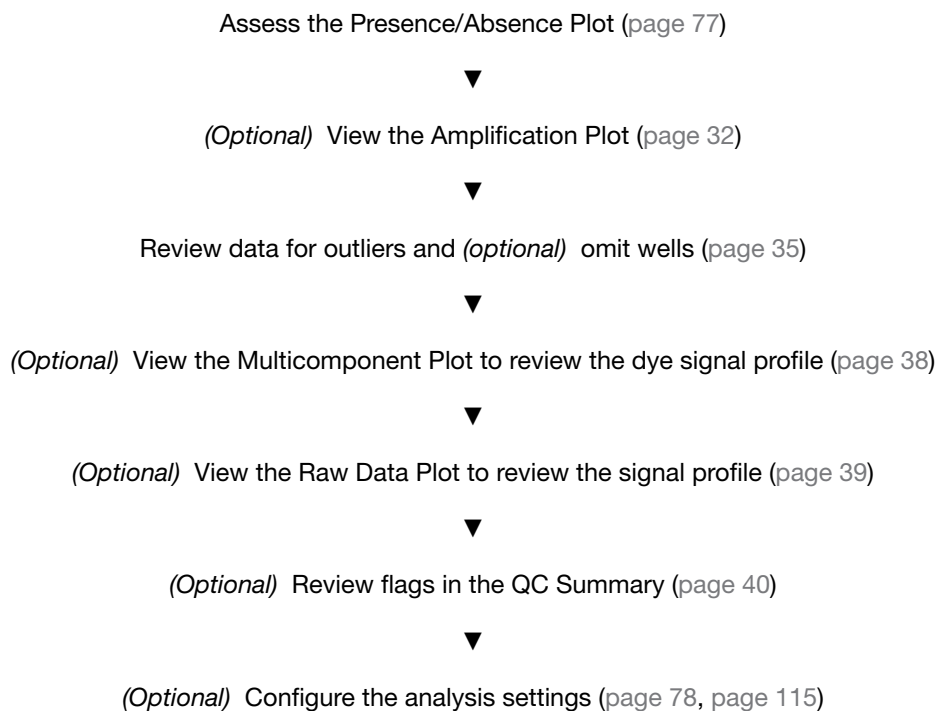

IMPORTANT! If you omit wells or configure the analysis settings, click **Analyze** to reanalyze the data.

Presence/Absence Plot overview

The Presence/Absence Plot displays the intensity of the fluorescence for each well.

Review the Presence/Absence Plot to confirm that amplification in the control wells is as expected and to review the calls for the unknown samples.

Table 22 Expected results for control reactions

Reaction type	Target	Result	Call
Negative control	IPC	Amplification	IPC Succeeded
	Target of interest ^[1]	No amplification	Negative control
NAC (blocked IPC)	IPC ^[2]	No amplification	Blocked IPC Control
	Target of interest	No amplification	Negative control

^[1] The target threshold is calculated from the negative control reactions.

^[2] The IPC threshold is calculated from the NAC reactions.

Table 23 Criteria for calls in unknown reactions

Target signal	IPC Signal	Call
Above the target threshold	Above or below the IPC threshold	Presence
Below the target threshold	Above the IPC threshold	Absence
Below the target threshold	Below the IPC threshold	Unconfirmed

View and assess the Presence/Absence Plot

You can view and assess the **Presence/Absence Plot** in the **Results** tab.

If no data are displayed in the **Results** tab, click **Analyze**.

1. In the **Results** tab, select **Presence/Absence Plot** from the dropdown list.
2. Click  to configure the plot, making the following selections:
 - **Target Reporter**: target defined for the target of interest
 - **Control Reporter**: target defined for the IPC
 - For the initial review of the **Presence/Absence Plot**, select the following options:
 - **Show Calls: All Calls**
 - **Show IPC**
 - **Show Controls**

The **Presence/Absence Plot** is displayed for data points selected in the plot settings. The data points for selected wells in the  **Plate Layout** or  **Well Table** are highlighted in the plot (see Figure 14).

3. Confirm that amplification in the negative and blocked IPC control wells is as expected. Use one of the following options:
 - Select control wells in the  **Plate Layout** or  **Well Table**, then confirm the location of the data points in the **Presence/Absence Plot**.
 - In the  **Well Table**, select **Group By ▶ Task**, then examine the wells with the Blocked IPC () and NTC () tasks. The C_q (C_t or C_{rt}) values should be Undetermined.
 - View the amplification plots for the negative controls (see Figure 15 and “Optimize display of negative controls in the Amplification Plot” on page 36).
4. In the **Presence/Absence Plot**, view the signal intensity and calls for the unknown samples. Use the plot settings (click ) to filter out the IPC results and control wells, or to select only one type of call.

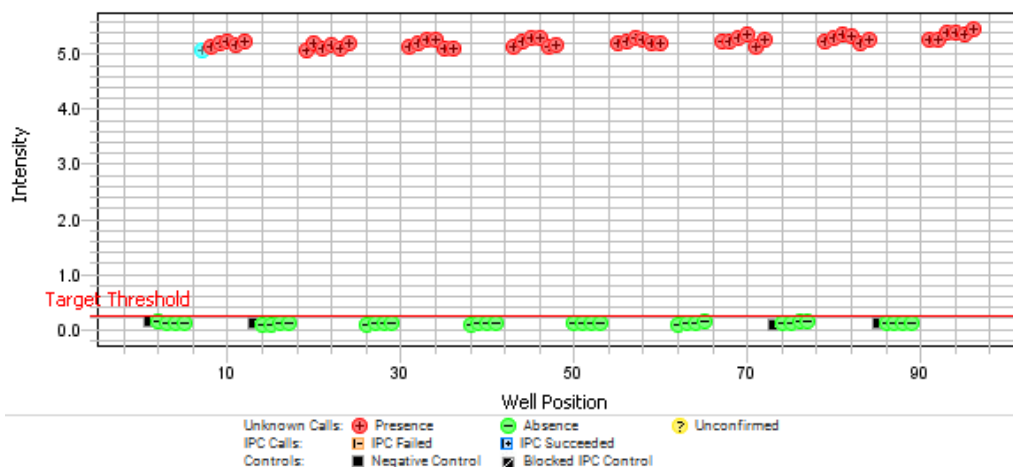


Figure 14 Example Presence/Absence Plot

The IPC results are not displayed in this example.

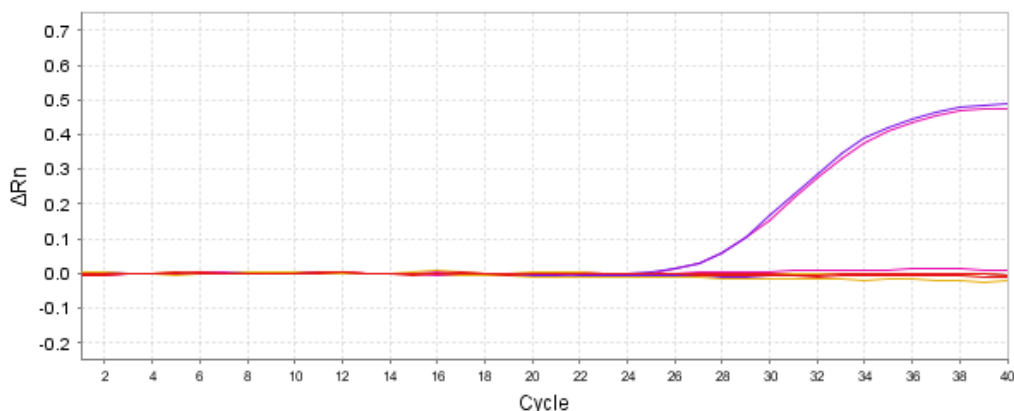


Figure 15 Example Amplification Plot for negative control and blocked IPC

Amplification of the IPC target (blue lines) is seen in the negative control wells but not the blocked IPC (NAC) wells. No amplification of the target of interest (red lines) is seen in either negative control or blocked IPC wells.

Call settings overview (presence / absence)

Use the **Call Settings** tab to edit the following settings:

- Data analysis settings
- Default call settings for assays without custom call settings
- Custom call settings for individual assays

Table 24 Options for data analysis settings (presence/absence experiments)

Data analysis setting	Description
Analyze Data from Post-PCR Read Only	Only post-PCR read data is used to determine calls.
Analyze Data from Pre-PCR Read and Post-PCR Read	The pre-PCR read is subtracted from the post-PCR read to determine calls.

Table 25 Options for call settings (presence/absence experiments)

Call setting	Description
Confidence Value	What confidence value is used to determine the target and IPC call thresholds. • A lower confidence value or more controls typically results in a lower calculated threshold. • A higher confidence value or fewer controls typically results in a higher calculated threshold.

For step-by-step instructions for adjusting the call settings, see the desktop software Help.



Set up, run, and review melt curve experiments

■ Melt curve experiments	80
■ Set up a melt curve experiment in the software	81
■ Set up melt curve reactions	82
■ Review results	82

Melt curve experiments

Overview

Use melt curve experiments to determine the melting temperature (T_m) of the amplification products of a PCR that used intercalating dyes.

Melting temperature (T_m) is the temperature at which 50% of the DNA is double-stranded and 50% is dissociated into single-stranded DNA. The melt curve of a single amplification product displays a single peak at the product's T_m . Multiple peaks in a melt curve experiment indicate additional amplification products, usually from non-specific amplification or formation of primer-dimers.

In the software, melt curve analysis is included in the default run method for any experiment type that uses intercalating dyes.

1. The software plots a melt curve based on the fluorescence of the dye with respect to change in temperature.
2. Using the melt curve, the software calculates the melting temperature (T_m).

Reaction types

Table 26 Reaction types for melt curve experiments

Reaction type (task)	Sample description
Unknown	Previously run PCR reactions that used intercalating dyes
No-template control (NTC/ Negative Control)	Previously run PCR reactions that used water or buffer Note: No DNA should be present in NTC wells.

Set up a melt curve experiment in the software

1. In the  **Home** screen, create or open a template.
 - In the  **New Experiment** pane, click **Create New Experiment** to create a new template.
 - In the  **Open Existing Experiment** pane, click **Open** to select and open an existing template.
2. In the **Properties** tab, enter the template information.
3. In the **Method** tab, adjust the reaction volume.
(Optional) Edit the ramp increment in the melt curve (see “Edit the ramp increment for the melt curve dissociation step” on page 106).
4. In the **Plate** tab (**Quick Setup**), assign plate attributes.
 - a. In the **Plate Attributes** pane, select the **Passive Reference** from the dropdown list.
5. In the **Plate** tab (**Quick Setup**), define and assign well attributes.
 - a. Select wells in the  **Plate Layout** or the  **Well Table**.
 - b. Assign samples and targets to selected wells.
 - Enter new sample and target names in the text fields.
 - Select previously defined samples and targets from the dropdown lists.

Note: New sample or target names entered in the **Quick Setup** subtab are automatically populated with default values for **Reporter** (FAM) and **Quencher** (NFQ-MGB) and assigned a task ( **Unknown**). Edit these values in the **Advanced Setup** subtab.

6. (Optional) In the **Plate** tab (**Advanced Setup**), assign tasks.
 - a. Select wells in the  **Plate Layout** or the  **Well Table**.
 - b. In the **Targets** table, select the checkbox of a target, then select a task from the dropdown list.


Reaction type	Task
Unknown (test sample)	
No-template control	

Set up melt curve reactions

Melt curve experiments are performed using previously amplified PCR products, usually at the end of the PCR run method. You can also use a plate from an intercalating dyes-based PCR run on another instrument.

Review results

Workflow: Review melt curve experiments

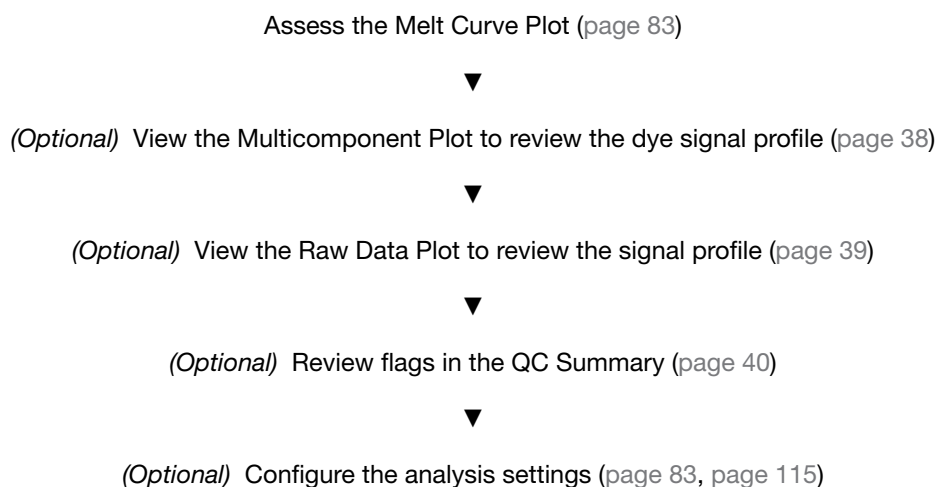

IMPORTANT! If you omit wells or configure the analysis settings, click **Analyze** to reanalyze the data.

Melt Curve Plot overview

The Melt Curve Plot displays the melt curve of the amplification products in the selected wells.

Review the Melt Curve Plot to confirm that the amplification products in a well display a single melting temperature (T_m). Multiple peaks in a melt curve indicate non-specific amplification or primer-dimer formation.

Table 27 Melt Curve plots

Plot	Description
Derivative Reporter vs. Temperature	Displays the derivative reporter signal in the y-axis as a function of temperature. The peaks in the plot indicate significant decrease in SYBR™ Green signal, and therefore the T_m of the amplified products. Use this plot to confirm a single T_m of the amplification products.
Normalized Reporter vs. Temperature	Displays the fluorescence signal from the reporter dye normalized to the fluorescence signal of the passive reference, as a function of temperature. You can use this plot to check the quality of the fluorescence data.

View and assess the Melt Curve Plot

You can view and assess the **Melt Curve Plot** in the **Results** tab.

If no data are displayed in the **Results** tab, click **Analyze**.

1. In the **Results** tab, select **Melt Curve Plot** from the dropdown list.
2. Click  to configure the plot, making the following selections:
 - **Plot Type:** **Derivative Reporter**
 - **Color:** **Sample**, **Target**, or **Well**
 - **Target:** **All** or a target of interest
 - *(For custom experiments with more than one Melt Curve stage)* Select the Melt Curve stage to view.

The **Melt Curve Plot** is displayed for the selected wells of the selected stage.

3. Review the plot for evidence of unexpected multiple peaks, which can indicate non-specific amplification or formation of primer-dimers.
4. Review the  **Well Table** for the calculated T_m in each well.

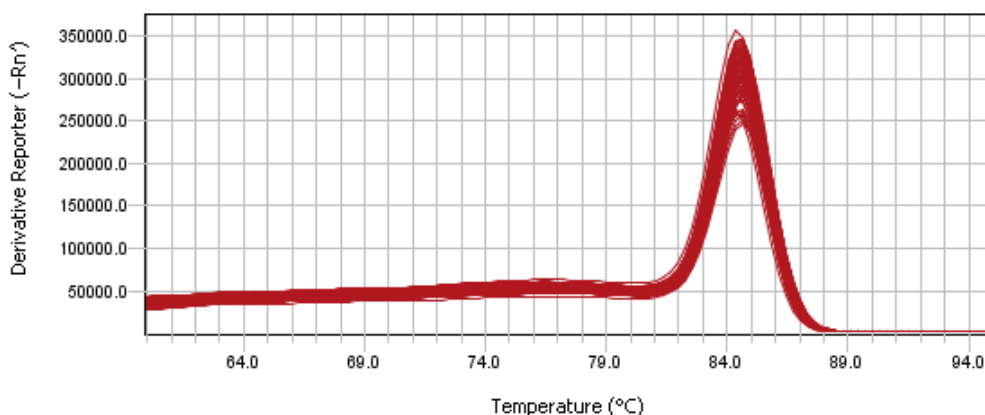


Figure 16 Example Melt Curve Plot

Melt curve settings overview

Use the **Melt Curve Settings** tab to enable or disable the Multi-Peak Calling function and adjust the detection levels for additional peaks, if needed.

- Enable or disable **Multi-Peak Calling**.

Multi-Peak Calling	Description
Enabled	– More than one PCR product is expected to amplify. – T_m will be determined for more than one peak.
Disabled	– A single PCR product is expected to amplify. – T_m will be determined for one peak.

- (For multi-peak calling only) Adjust the detection levels for additional peaks.

Option	Description
Peak level relative to the dominant peak (%)	Specify a fractional-level value as the additional peak detection threshold. The detected peaks are measured relative to the height of the tallest peak, which has a perfect fractional level of 100%. The default value is 10%. For example, set a fractional-level detection threshold value at 40, then peaks above 40% of the tallest peak are reported, and peaks below 40% are regarded as noise.
Peak Calling Threshold	Specify an absolute fluorescence-level value as the peak calling threshold. The absolute fluorescence is measured on the derivative reporter (-dRn') axis. Only peaks that appear above the peak calling threshold will be detected. For example, set a fluorescence-level value at 90,000, then peaks with fluorescence above 90,000 are reported, and peaks below 90,000 are regarded as noise.

For step-by-step instructions for adjusting the melt curve settings, see the desktop software Help.



Use the software with the Security, Auditing, and E-signature (SAE) Administrator Console

Overview of the Security, Auditing, and E-signature Administrator Console

The software is compatible with the Security, Auditing, and E-signature Administrator Console.

When you install the QuantStudio™ Design and Analysis Desktop Software, you have the option to install the SAE Administrator Console.

The information in this chapter is for the security, audit, and e-signature tasks that are performed in the QuantStudio™ Design and Analysis Desktop Software.

For information about the Security, Auditing, and E-signature Administrator Console and the tasks that are performed in the console, see *SAE Administrator Console v2.1 User Guide* (Pub. No. MAN0028414).

Table 28 Functions of the Security, Auditing, and E-signature Administrator Console

Feature	Description
System security	Controls user access to an application. A default user account assigned the Administrator role is provided at installation. You can set up additional SAE user accounts and permissions.
Auditing	Tracks actions performed by users and changes to the SAE settings. Some actions are automatically audited silently. The following audit functions can be set up: • Audit changes to specific objects and specify the audit mode. • Generate reports for audited user actions and SAE function changes. • Generate reports for software or instrument actions and runs.
Electronic signature (e-signature)	Determines if users are required to fulfill signature requirements before performing specific functions. You can perform the following functions: • Configure e-signature so that a user can perform an action only if the associated data are signed. • Configure each e-signature event to require one or multiple signatures and to require users with specific roles to sign.

Enable the system security function

1. In the QuantStudio™ Design and Analysis Desktop Software, click **Tools ▸ Security ▸ Settings**.
2. In the **Login** dialog box, enter your user name and password, then click **Settings**.
3. In the **SAE Settings** dialog box, select the connection option.

Option	Description
Local server radio button	Use this option when the SAE Administrator Console and the QuantStudio™ Design and Analysis Desktop Software are installed on the same computer.
Remote server radio button	Use this option when the SAE Administrator Console and the QuantStudio™ Design and Analysis Desktop Software are installed on separate computers.

4. If a remote server was selected, enter the IP address of the computer that is running the SAE Administrator Console in the **Server IP address** field.
Port number 8443 must be used. Do not edit this value. Ensure that this port is open on the network.
Depending on the way that the network is configured, the IP address of the computer might change. This affects the connection between the SAE Administrator Console and the QuantStudio™ Design and Analysis Desktop Software. A DHCP network connection with a reserved IP address is recommended. A static IP address can also be used.
5. Click **Log In**.
6. In the **SAE Settings** dialog box, select the **Enable SAE settings** checkbox, then click **Save**.
7. Click **OK**.

Disable the system security function

1. In the QuantStudio™ Design and Analysis Desktop Software, click **Tools ▸ Security ▸ Settings**.
2. In the **SAE Settings** dialog box, deselect the **Enable SAE settings** checkbox, then click **Save**.
3. Click **Yes** to confirm.
4. Click **OK** to confirm.

Files that were created or edited with the security function enabled can be opened and edited after the security function is disabled. The audit and e-signature information for these files is not available after the security function is disabled. Any edits that are made to the files do not have any audit information collected. E-signatures are not collected.

The audit and e-signature information that was collected before the security function was disabled is retained. It can be accessed if the security function is enabled again.

Overview of accounts

You must sign in to the software when security is enabled. A user name and a password are required.

Accounts are set up and edited in the Security, Auditing, and E-signature Administrator Console.

Accounts are associated with permissions for specific tasks in the software.

For more information about setting up and editing accounts, see *SAE Administrator Console v2.1 User Guide* (Pub. No. MAN0028414).

In the QuantStudio™ Design and Analysis Desktop Software, individual users can change their own password. Administrators can change passwords for any account in the Security, Auditing, and E-signature Administrator Console.

Depending on how the account is set up in the Security, Auditing, and E-signature Administrator Console, an individual user might need to reset their password the first time they log in to the QuantStudio™ Design and Analysis Desktop Software.

Log in to QuantStudio™ Design and Analysis Desktop Software

If security is enabled in the QuantStudio™ Design and Analysis Desktop Software, the **Login** dialog box is displayed when the software is launched.

Ensure that the administrator has set up an account in the Security, Auditing, and E-signature Administrator Console.

1. In the **Login** dialog box, enter your user name in the **User Name** field, then enter your password in the **Password** field.
2. Click **Log In**.

Depending on how the account is set up in the Security, Auditing, and E-signature Administrator Console, an individual user might need to reset their password the first time they log in to the QuantStudio™ Design and Analysis Desktop Software.

Change your password in QuantStudio™ Design and Analysis Desktop Software

1. In the QuantStudio™ Design and Analysis Desktop Software, click **Tools** ▶ **Security** ▶ **Change Password**.
2. Enter your password in the **Old Password** field.
3. Enter your new password in the **New Password** field and the **Confirm Password** field.
4. Click **OK**.

Log out of QuantStudio™ Design and Analysis Desktop Software

In the QuantStudio™ Design and Analysis Desktop Software, click **Tools ▶ Security ▶ Log Out**.

Default roles and permissions

IMPORTANT! Permissions apply to all of the user accounts that are assigned to the role.

Table 29 Default roles and permissions

Permissions		Roles		
Category	Function	Administrator	Scientist	Technician
Setup	Create new template	Yes	Yes	Yes
	Edit template properties	Yes	Yes	Yes
	Edit run method	Yes	Yes	No
	Edit analysis settings	Yes	Yes	No
	Omit wells for analysis	Yes	Yes	No
	Define new target/assay	Yes	Yes	No
	Define new sample	Yes	Yes	Yes
	Edit sample	Yes	Yes	Yes
	Delete sample	Yes	Yes	Yes
	Assign targets/assays	Yes	Yes	Yes
	Assign samples	Yes	Yes	Yes
	Edit passive reference	Yes	Yes	Yes
Run	Start a run from the desktop software	Yes	Yes	Yes
	Stop a run from the desktop software	Yes	Yes	Yes
Target library	Create targets	Yes	Yes	Yes
	Edit targets	Yes	Yes	Yes
	Delete targets	Yes	Yes	No
Analysis settings library	Create analysis settings	Yes	Yes	Yes
	Edit analysis settings	Yes	Yes	Yes
	Delete analysis settings	Yes	Yes	No

Table 29 Default roles and permissions (continued)

Permissions		Roles		
Category	Function	Administrator	Scientist	Technician
Dye library	Create a custom dye	Yes	Yes	Yes
	Delete a custom dye	Yes	Yes	No
Preferences	Edit system preferences	Yes	Yes	No
Edit file save destination	Edit save as destination	Yes	Yes	No
	Edit export destination	Yes	Yes	No
	Edit report destination	Yes	Yes	No
Security configuration	Log into timed-out user sessions	Yes	No	No
	Perform e-signing	Yes	Yes	Yes
Security configuration (Security, Auditing, and E-signature Administrator Console)	Configure security and auditing	Yes	No	No
	View action records	Yes	No	No
	View system configuration records	Yes	No	No
	View application object records	Yes	No	No
	View instrument run records	Yes	No	No

View and manage the audit records

View the object audit records

The object audit records are maintained within the template file or the data file. The audit records are not displayed in the SAE Administrator Console.

1. In the QuantStudio™ Design and Analysis Desktop Software, open a template file or a data file.
2. Select the **Audit** tab.
The audit records for the file are displayed. There are two panes, a **Summary** pane and a **Details** pane. The **Summary** pane lists the audit records. Each audit record corresponds to an update of the file. The **Details** pane lists each individual edit to a file.
3. In the **Summary** pane, select a row that corresponds to an audit record.
The details for the specific audit record are displayed in the **Details** pane.
4. (Optional) Click  **Action ▶ Show Filter** to display the filter pane at the top of the screen.
If the filter pane is displayed, the button is **Hide Filter**.

5. Select the **Filter by** checkbox to filter the audit records.

The audit records can be filtered by the following items:

- Date range
- User name
- Full name
- Record type
- Record name
- Reason
- Action

IMPORTANT! When entering a value in a field, the text must match the exact text in the audit record. It is case-sensitive and the punctuation must match. If the text does not match the audit record, no results are displayed.

6. Click **Refresh** to update the results.

The results display the audit records that correspond to the filters. The audit records that correspond to the filters are displayed in the **Summary** pane. Select a row to view the details in the **Details** pane.

Some of the filters are based on the details. For example, if you filter by an action, all of the audit records that contain this action in the details are displayed in the **Summary** pane.

Generate an object audit report

All of the audit records are included in the report. This applies even if filters were applied in the software.

1. In the QuantStudio™ Design and Analysis Desktop Software, open a template file or a data file.
2. Select the **Audit** tab.
3. Click  **Action ▶ View Report**.
The **Print Preview** dialog box is opened and it displays a preview of the audit report.
4. Click **Save As PDF**, select a folder location, enter a file name, then click **Save**.
5. Click **Print Report** to print the file.

Export the audit summary

1. In the QuantStudio™ Design and Analysis Desktop Software, open a template file or a data file.
2. Select the **Audit** tab.
3. In the **Summary** tab, select a row that corresponds to an audit record.
4. Click  **Action ▶ Export Summary**.
5. In the **Save** dialog box, select a folder location, enter a file name, then click **Save**.

The summary is exported in a CSV file format.

Export the audit details

1. In the QuantStudio™ Design and Analysis Desktop Software, open a template file or a data file.
2. Select the **Audit** tab.
3. In the **Summary** tab, select a row that corresponds to an audit record.
4. Click  **Action ▶ Export Details**.
5. In the **Save** dialog box, select a folder location, enter a file name, then click **Save**.

The details are exported in a CSV file format. The information that is included in the file is from the **Details** pane that is associated with the selected audit record.

View and manage the e-signatures

Sign the data

1. In the QuantStudio™ Design and Analysis Desktop Software, open a template file or a data file.
2. Click **Tools ▶ Security ▶ Sign data**.
3. In the **Electronic Signature** dialog box, select a meaning from the **Meaning** dropdown list.
The meanings that are available in the dropdown lists correspond to the meanings that were set up in the Security, Auditing, and E-signature Administrator Console.
4. Enter your user name and password, then click **Sign Data**.

The electronic signature is available in the **e-sig** tab.

View the e-signatures

1. In the QuantStudio™ Design and Analysis Desktop Software, open a template file or a data file.
2. Select the **e-Sig** tab.
All of the e-signatures for the file are displayed.
3. (Optional) Select a row to view a preview of the experiment report in the bottom pane.
4. (Optional) Save or print the experiment report.
 - Click **Save As PDF**, select a folder location, enter a file name, then click **Save**.
 - Click **Print Report**, then select the print settings.
5. (Optional) Click  **Action ▶ Print e-Signatures**.

6. (Optional) In the **Print Preview** dialog box, save or print the e-signature history report.

- Click **Save As PDF**, select a folder location, enter a file name, then click **Save**.
- Click **Print Report**, then select the print settings.



Instrument overview


■ Power on the instrument	93
■ View run history and delete or transfer files from the instrument	93
■ Load and unload the plate in the instrument	94

Power on the instrument

1. Touch anywhere on the touchscreen to determine if the instrument is in sleep mode. If the home screen is displayed, the instrument is already powered on.
2. If the home screen does not display, power on the instrument by pressing the switch on the rear panel.

If left unattended (for about two hours), the instrument automatically enters sleep mode (enabled by default) to conserve power.

Note: To customize the sleep mode setting, touch  **Settings ▶ Instrument Settings ▶ Sleep Mode.**

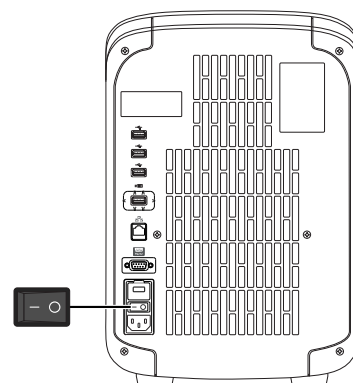

View run history and delete or transfer files from the instrument

In the home screen, touch  **Settings ▶ Run History.**

- Touch an individual run record to view its details, then complete one of the following actions:
 - Touch **Delete** to delete the run record.
 - Touch **Transfer** to export the run data.
- Touch **Manage** to select multiple run records for simultaneous viewing, deletion, or transfer.

Note:

- Guests (users not signed-in) can only view guest run records.
- Users signed into their instrument profiles can also view their own run records.
- Administrators can view all run records.

Note: If the connection between the instrument and the desktop software is interrupted during the run, the instrument still completes the run. However, the run data (EDS file) must be transferred from the instrument to the desktop software using a USB drive or a network drive.


Load and unload the plate in the instrument

Load and unload a plate in the QuantStudio™ 1 Real-Time PCR Instrument



CAUTION! Use flat caps for 0.2-mL tubes. Rounded caps can damage the heated cover.

1. Load the plate.
 - a. Open the instrument drawer.
 - b. Load the plate onto the plate adapter so that the following criteria are met.
 - Well A1 of the plate is in the top-left corner of the plate adapter.
 - The barcode faces the front of the instrument.

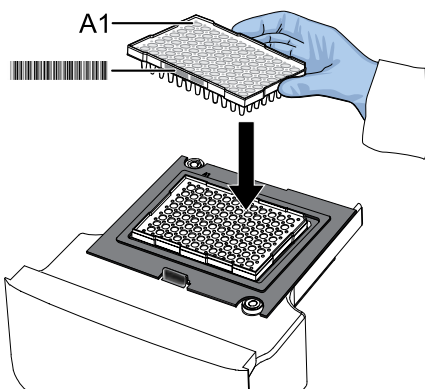

Note: Do not remove the black plate adapter before loading a plate or tube strips. If used, tube strips can fit loosely in the adapter, but the heated cover will apply the appropriate pressure to seat the tube strips securely in the adapter.

- c. Close the instrument drawer.
2. When the run ends, unload the plate.
 - a. Open the instrument drawer.
 - b. Remove the plate.

- c. Close the instrument drawer.




CAUTION! PHYSICAL INJURY HAZARD. During instrument operation, the plate temperature can reach 100°C. Allow it to cool to room temperature before handling.

Load and unload a plate in the QuantStudio™ 3 Real-Time PCR Instrument or QuantStudio™ 5 Real-Time PCR Instrument



CAUTION! Use flat caps for 0.2-mL tubes and 0.1-mL tubes. Rounded caps can damage the heated cover.

1. Load the plate.
 - a. Touch to eject the instrument drawer.
 - b. Load the plate onto the plate adaptor so that the following criteria are met.
 - Position well A1 of the plate in the top-left corner of the plate adapter.
 - Ensure that the barcode faces the front of the instrument.

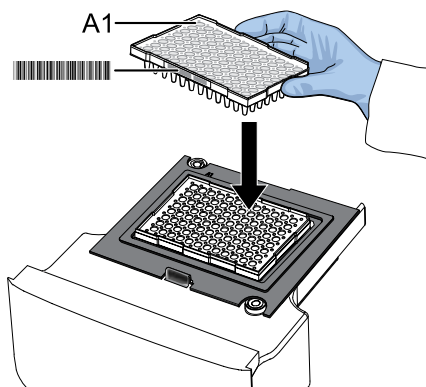

IMPORTANT! The instrument should be used by trained operators who have been warned of the moving parts hazard.

Note: (For 96-well 0.2-mL blocks only) Do not remove the black plate adapter before loading a plate or tube strips. If used, tube strips can fit loosely in the adapter, but the heated cover will apply the appropriate pressure to seat the tube strips securely in the adapter.

Note: The 384-well and 96-well Fast (0.1-mL) block configurations do not require a plate adapter.

- c. Touch to close the instrument drawer.
2. When the run ends, unload the plate.
 - a. Touch to eject the instrument drawer.
 - b. Remove the plate.

- c. Touch  to close the instrument drawer.



CAUTION! PHYSICAL INJURY HAZARD. During instrument operation, the plate temperature can reach 100°C. Allow it to cool to room temperature before handling.

Note: If the instrument does not eject the plate, contact Support.



Alternative procedures to set up a template

■ Set up a custom experiment	97
■ Assign samples using a sample definition file	99
■ Assign samples and targets using plate setup files	100
■ Assign targets, samples, and biological replicate groups from an XLS file	101
■ Create new EDT files using existing EDT and EDS files	102

Set up a custom experiment

Custom experiments overview

Custom experiment setup is required for assays that use multiple PCR stages, such as TaqMan™ Mutation Detection Assays. A custom experiment also allows flexibility for secondary analysis.

The default settings for custom experiments are that of a standard curve experiment, but most settings are editable.

Table 30 Default settings in custom experiments

Setting	Default
Run method (thermal protocol)	Equivalent to standard curve experiment default
Tasks	Unknown  Negative control  Standard 
C _t settings	Baseline threshold equivalent to standard curve experiment default
Flag settings	QC flags on No automatic omissions
Auto Export	Off

Table 31 Editing options in custom experiments

Setting	Description
Run method (thermal protocol)	- Ramp rates can be edited within software limits. - Data collection can be enabled at any step and during any ramp within a melt stage. - Multiple instances of any type of stage can be added, with exceptions noted. - Any stage can be added at any point in a run method, with exceptions noted. Noted exceptions: - Only one infinite hold (must be added at the end). - Only one pre-PCR read and one post-PCR read stage. If both exist in a run method, the pre-read must be before the post-read. For example, the following order is valid: melt-PCR-Pre-Read-Melt-PCR.
Analysis settings	Editable
Flags	Editable
Optical filters	Editable

Set up a custom experiment in the software

- In the  **Home** screen, create or open a template.
 - In the  **New Experiment** pane, click **Create New Experiment** to create a new template.
 - In the  **Open Existing Experiment** pane, click **Open** to select and open an existing template.
- In the **Properties** tab, enter the template information.
- In the **Method** tab, edit the default run method according to the experiment requirements.
- In the **Plate** tab (**Quick Setup**), assign plate attributes.
 - In the **Plate Attributes** pane, select the **Passive Reference** from the dropdown list.
- In the **Plate** tab (**Quick Setup**), define and assign well attributes.
 - Select wells in the  **Plate Layout** or the  **Well Table**.
 - Assign samples and targets to selected wells.
 - Enter new sample and target names in the text fields.
 - Select previously defined samples and targets from the dropdown lists.

Note: New sample or target names entered in the **Quick Setup** subtab are automatically populated with default values for **Reporter** (FAM) and **Quencher** (NFQ-MGB) and assigned a task ( **Unknown**). Edit these values in the **Advanced Setup** subtab.

6. (Optional) In the **Plate** tab (**Advanced Setup**), assign tasks.
 - a. Select wells in the  **Plate Layout** or the  **Well Table**.
 - b. In the **Targets** table, select the checkbox of a target, then select a task from the dropdown list.
7. (Optional) In the **Plate** tab (**Advanced Setup**), define and assign biological replicate groups (see “Define and assign biological replicate groups” on page 110).

Assign samples using a sample definition file

About sample definition files

Import sample information from a sample definition file to include in the plate setup. A sample definition file is a CSV file or a TXT file that contains the following setup information:

- Sample name
- (Optional) Custom sample properties

Create a sample definition file

1. In a spreadsheet program, create the following column headers:
 - Well
 - Sample Name
 - (Optional) Column header names for up to 32 user-defined custom fields (for example, **Custom 1**, **Custom 2**, etc.)
2. Enter the well number and sample name in the appropriate columns.
3. (Optional) Enter the custom properties for the sample.
4. Save the file as a tab-delimited text file (TXT) or a comma-separated values file (CSV).

Import sample information from a sample definition file

Example setup files are provided with the software in:

`<drive> : \Program Files (x86)\Applied Biosystems\QuantStudio Design and Analysis Software\examples\User Sample Files,`

where `<drive>` is the drive on which the software is installed.

1. In an open experiment, select **File ▶ Import Plate Setup**.
2. Click **Browse**, navigate to a sample definition text file, then click **Select**.
3. Click **Apply**.
4. If the experiment already contains plate setup information, the software prompts for the replacement of the plate setup with the data from the file. Click **Yes** to replace the plate setup information.

The samples appear in the **Samples** table for the experiment. All samples and well assignments in the experiment are replaced with those in the file. If defined, the custom sample properties also appear in the ≡ **Well Table** of the **Results** tab and in the ≡ **Plate Layout** tooltips in both the **Plate** and **Results** tabs. The custom fields can be exported with the results data.

Note: To modify custom sample properties information, edit the custom fields in the sample definition file and import the file again. All of the sample information in the experiment is replaced with the information in the new file.

Assign samples and targets using plate setup files

About plate setup files

Plate setup files contain setup information such as the well number, sample name, sample color, target name, dyes, and other reaction plate contents.

Plate setup files can be exported from previously run experiments. For instructions on exporting an experiment, see “Export experiments or results” on page 42.

Import plate setup data

Import the plate setup for a new experiment from an exported file with one of the following formats:

- EDS – EDS file format
- EDT – user-created and system templates files format
- TXT – text format
- XML – XML format
- CSV – comma separated values format
- SDT – Sequence Detection System (SDS) template files format
- SDS – 7900 v2.4 format

Note: Import plate setup information from a 96-well plate into a 384-well plate, provided that the sample file is a TXT file.

IMPORTANT! The file must contain only plate setup data and it must match the experiment type.

1. In the **Plate** tab, select **File ► Import Plate Setup**.
2. Click **Browse**, navigate to and select the file to import, then click **Select**.

Example setup files are provided with the software in:

<drive> : \Program Files (x86)\Applied Biosystems\QuantStudio Design and Analysis Software\examples\User Sample Files,
where <drive> is the drive on which the software is installed.

3. Click **Apply**.
The setup data from the selected file is imported into the open experiment.

Assign targets, samples, and biological replicate groups from an XLS file

For wells with single targets, you can paste assignment information from an XLS file into the plate layout of the desktop software.

An example copy and paste file is provided with the software in:

<drive> : \Program Files (x86)\Applied Biosystems\QuantStudio Design and Analysis Software\examples\User Sample Files.

where <drive> is the drive on which the software is installed.

1. In the custom properties tab of the example Microsoft™ Excel™ file, ensure that the **Well** column is sorted in order 1 through 96, then select the **Well** column and the **Sample Name** column, *including* headers.
2. In the **Plate** tab of the software, click ≡ **Well Table**, then ensure that the well numbers are in order from 1 through 96.
3. In the ≡ **Well Table**, hover the mouse in the first cell underneath the **Sample** header (adjacent to A1), right-click, then select either **Paste** or **Paste only samples**.
Any of the columns not copied are treated as NULL values for those columns.

Create new EDT files using existing EDT and EDS files

About experiment templates

Use templates to create experiments with the same parameters or with pre-existing settings. Experiments can be saved as unlocked or locked (password-protected) templates.

You can save the following information in an experiment template (EDT) file:

- Plate setup information (defined sample and targets or SNP assays, plate assignment of samples and targets or SNP assays)
- Reagent information
- Run method (thermal protocol)
- Analysis settings

Example templates are provided with the software in:

```
<drive>:\Program Files (x86)\Applied Biosystems\QuantStudio Design & Analysis  
Software\templates,
```

where <drive> is the drive on which the software is installed.



Detailed procedures to create or edit a method

■ Adjust method parameters	103
■ Set up advanced temperature zones (Auto Delta and VeriFlex™ Zones)	104
■ Add or adjust a pause step	105
■ Select optical filters	105
■ Edit the ramp increment for the melt curve dissociation step	106

Adjust method parameters

For an overview of the method as it is graphically represented, see “Method elements” on page 19.

- In the **Method** tab, click a method parameter field to edit the following information:
 - Reaction volume
 - Temperature ramp rate
 - Step temperature
 - Step hold time
 - Number of cycles
- Click-drag  to increase or decrease a step temperature.
- Click  to switch data collection on or off at each step.
Data Collection On enables analysis of data that is collected throughout the PCR, for real-time analysis and troubleshooting.
- (QuantStudio™ 3 Real-Time PCR Instrument and QuantStudio™ 5 Real-Time PCR Instrument only)
Click  to configure settings for **Auto Delta** or VeriFlex™ Zones for individual steps (see “Set up advanced temperature zones (Auto Delta and VeriFlex™ Zones)” on page 104).

Note: In melt curve stages, **Advanced Settings** are not applicable.

An **A** or **V** is displayed in the PCR stage when **Auto Delta** or VeriFlex™, respectively, is enabled.

- Click  to configure pause settings.
- Adjust the heated cover temperature via the instrument settings (see the instrument user guide).



- See the Help to learn more about adjusting the following parameters:
 - Adding or subtracting a stage
 - Adding or subtracting a step from a stage
 - Configuring optical filter settings

Set up advanced temperature zones (Auto Delta and VeriFlex™ Zones)

(QuantStudio™ 3 Real-Time PCR System and QuantStudio™ 5 Real-Time PCR System only) Configure settings for **Auto Delta** and **VeriFlex™ Zones**.

- **Auto Delta**—Incremental increase or decrease of a cycle's temperature or hold time for a step in a cycling stage (not applicable for **Hold** or **Infinite Hold** stages).
- **VeriFlex™ Zones**—Independent temperature zones within 5°C of adjacent zones.
 - QuantStudio™ 3 Real-Time PCR Instrument: 3 zones
 - QuantStudio™ 5 Real-Time PCR Instrument: 6 zones

Note: VeriFlex™ Zones temperature settings are not available for 384-well blocks.

1. In the **Method** tab, click  **Advanced Settings** in a step.

Note: Any changes apply only to the step in which you clicked.

2. Configure either the **VeriFlex™ Zones** or **Auto Delta** for the selected step.

- Select **VeriFlex™**, then enter a temperature for each zone.

Note: In the **Plate** tab, the VeriFlex™ Zones display on the plate layout.

- To view setting details, hover over the **V** in each zone.
 - To hide the display of zones, click  **Action ▶ Hide VeriFlex™ Zones**.
-

- Select **Auto Delta**, then enter a starting cycle, temperature, and time.

3. Click **Save**.



Add or adjust a pause step

1. In the **Method** tab, click  in the step.
2. Select **Pause**.
3. Enter the cycle after which the pause should occur.
4. Enter a pause temperature between 4°C and 99.9°C.



CAUTION! PHYSICAL INJURY HAZARD. During instrument operation, the plate temperature can reach 100°C. If you want to access the plate during a run pause, enter room temperature as the pause temperature and allow the plate to cool to room temperature before handling.

5. Click outside of the pause dialog box to return to the method.
6. (Optional) To remove a pause, click , then deselect **Pause**.

Select optical filters

The need to edit optical filter settings is rare, and it is intended for advanced or custom uses only. For information on the dyes read by each filter, see the instrument user guide.

Use the optical filters settings feature to select a filter set to match the profile of a custom dye.

1. In the **Method** tab, select  **Action ▶ Optical filter settings**.
The excitation (x) and emission (m) wavelengths that correspond to each filter are shown on the screen.
2. Select the check boxes to enable or disable filters.
A **Melt Curve Filter** table is accessible if the method contains a melt curve stage. Otherwise, use the **PCR Filter** table to select optical filters.
3. (Optional) Click **Revert to Defaults** to reset filters.
4. Click **Close**.



Edit the ramp increment for the melt curve dissociation step

In the **Method** tab, you can perform the following tasks to edit the ramp increment for the melt cure dissociation step.

- Select the ramp increment method for the dissociation step (located under the graphical representation of the thermal protocol).

Option	Description
Continuous (default)	Continuously increases the temperature by the ramp increment (°C/sec).
Step and Hold	Increases the temperature by the ramp increment (°C), then holds at that temperature for the specified time.
No. of Data Points per Degree	Increases the temperature by the ramp increment (°C) and collects the specified number of data points per degree increased.

- (For all options) Edit the temperature ramp increment.
 - a. Click the ramp increment element in the **Dissociation** step.
 - b. Enter a value or use the up/down arrows (default is 0.15°C/s).
- (Step and Hold only) Edit the hold time after each temperature increase.
 - a. Click the time field next to **Step and Hold**.
 - b. Enter a value or use the up/down arrows (default is 5 seconds).
- (No. of Data Points per Degree only) Edit the number of data points to be collected with each degree increase.
 - a. Click the number of data points element in the **Dissociation** step.
 - b. Enter a value or use the up/down arrows (default is 10 data points).



Detailed procedures to set up plate / well details and libraries

■ Assign well attributes (Quick Setup subtab)	107
■ Define and assign well attributes (Advanced Setup subtab)	107
■ Assign a task to wells	108
■ Define and set up standard dilutions	109
■ Define and assign biological replicate groups	110
■ Sample, target, and SNP assay libraries	111

Assign well attributes (Quick Setup subtab)

In the **Quick Setup** subtab, assign well attributes by direct entry into text fields or by selecting user-defined samples and targets or SNP assays from dropdown lists.

1. In the **Plate** tab, click **Quick Setup**.
2. Select plate wells in the  **Plate Layout** or the  **Well Table**.
3. Assign the well attributes for the selected wells.
 - Into the text fields, enter the names of the new sample and the new target or SNP assay.
 - From the dropdown lists, select a user-defined sample and target or SNP assay.
For more information about defining or importing samples and targets or SNP assays, see “Define and assign well attributes (Advanced Setup subtab)”.

Note: In the **Advanced Setup** subtab, change the default selections for the reporter and quencher dyes and for tasks (see “Assign a task to wells” on page 108).

4. (Optional) Enter comments for the selected wells.

Define and assign well attributes (Advanced Setup subtab)

In the **Advanced Setup** subtab, define or import samples and targets or SNP assays, then assign well attributes.

1. In the **Plate** tab, click **Advanced Setup**.
2. In the **Samples** table, define samples (see “Define samples in the Samples table” on page 112).

3. In the **Targets** or **SNP Assays** table, define targets or SNP assays, then select detection tasks.
 - a. Define targets or SNP assays (see “Define targets in the Targets table” on page 113 or “Define SNP assays in the SNP Assays table” on page 114, respectively).
 - b. Select a detection task from the **Task** column dropdown list (see “Assign a task to wells” on page 108).
4. Assign well attributes.
 - a. Select plate wells in the  **Plate Layout** or the  **Well Table** (see “Select plate wells” on page 21).
 - b. Select the checkbox of a defined sample.
 - c. Select the checkbox of a defined target or SNP assay.

Assign a task to wells

1. In the **Plate** tab, click **Advanced Setup**.
2. Select plate wells in the  **Plate Layout** or the  **Well Table** (see “Select plate wells” on page 21).
3. In the **Targets** or **SNPs** table, select the check box of a target or SNP assay.
4. Select a detection task from the **Task** column dropdown list.

Detection tasks for targets and SNP assays

Task	Description
 Unknown (default)	The well contains test samples with unknown genotype.
 Negative Control / No template control	The well contains water or buffer instead of sample.
 Standard ^[1]	The well contains samples with known standard quantities. Note: For a standard detection task, enter the standard quantity in the quantity column.
 Positive Control Allele 1 / Allele 1 ^[2]	The well contains samples homozygous for allele 1.
 Positive Control Allele 2 / Allele 2 ^[2]	The well contains samples homozygous for allele 2.
 Positive Control Allele 1 / Allele 2 ^[2]	The well contains samples heterozygous for allele 1 and 2.

(continued)

Task	Description
 Internal positive control ^[3]	The PCR reaction contains a short synthetic DNA template to distinguish between true negative results (that is, the target is absent in the samples) and negative results caused by PCR inhibitors, incorrect assay setup, or reagent or instrument failure.
 Blocked IPC ^[3]	The well contains an IPC blocking agent, which blocks amplification of the IPC.
NAC – No amplification control ^[3]	The PCR reaction contains an IPC blocking agent instead of sample. No amplification should occur in negative control-blocked IPC wells because the reaction contains no sample and amplification of the IPC is blocked.

^[1] For standard curve and relative standard curve experiments only.

^[2] For genotyping experiments only.

^[3] For presence/absence experiments only.

Define and set up standard dilutions

Note: This information is applicable for standard curve and relative standard curve experiments only.

1. In the **Plate** tab, select  **Action ▶ Define and Set Up Standards.**
2. Select **Singleplex** or **Multiplex** from the **Model** dropdown list.
3. (Optional) Select the target from the dropdown list.
4. Enter the parameters for the dilution series.
 - **Number of dilution points**—5 recommended
 - **Number of replicates**—3 recommended
 - **Starting Quantity**—The highest or lowest standard quantity, without units.

Note: The quantity can be expressed as copies, copies/μL, ng/μL, or as relative dilutions.

Note: Use **E** to indicate the exponent number using scientific notation. For example, to indicate 1.23×10^4 , enter 1.23E4.

- **Serial Factor**

Note: The serial factor calculates quantities for all standard curve points.

- Starting quantity is the highest value—Select 1:10 to 1:2.
- Starting quantity is the lowest value—Select 2× to 10×.

Note: The **Standard Curve Preview** y-axis values are calculated from the starting quantity and serial factor. Actual results may differ from the preview.

5. Select and arrange the wells to use for the standards.
 - Select **Automatically Select Wells for Me**.
 - Select **Let Me Select Wells**, then select wells using the displayed plate layout.
6. Select to arrange the standards in **Columns** or **Rows**.
7. (Optional) Click **Reset** to revert to default values.
8. Click **Apply**, then click **Close** to return to the **Plate** tab.

Assign the standard dilutions manually

Note: Applicable for standard curve and relative standard curve experiments only.

1. In the **Plate** tab, select wells in the  **Plate Layout** or  **Well Table**.
2. Select the check box for the target, select  from the **Task** dropdown list, then enter a quantity.
3. Repeat to complete the standard dilution series.

Define and assign biological replicate groups

Biological Replicate Groups can be used in standard curve, relative standard curve, comparative C_t , and custom experiments.

1. In the **Plate** tab, click **Advanced Setup**.
2. Define **Biological Replicate Groups**.
 - a. In the **Biological Replicate Groups** table, click  **Add**.
 - b. (Optional) Click a cell to edit color, name, or comments.
 - c. (Optional) Click  to delete a biological replicate group from the table.
3. Assign **Biological Replicate Groups**.
 - a. Select plate wells in the  **Plate Layout** or the  **Well Table** (see “Select plate wells” on page 21).
 - b. In the **Biological Replicate Groups** table, select the check box of a biological replicate group.



Sample, target, and SNP assay libraries

Libraries overview

Libraries contain saved information to reuse in future templates.

The following libraries are available in the software:

- Dye Library
- Sample Library
- Target Library
- SNP Assay Library
- Analysis Settings Library

To access the libraries, use the **Tools** menu or the **Plate** tab.

- In the menu bar, select **Tools ▶ (Library of choice)**.
- In the **Plate** tab of an open EDT or EDS file, click **Advanced Setup**, then select  **Action ▶ Import from Library**.

Permissions for the libraries

If security is enabled in the software, the following permissions are applicable:

- Targets library
 - **Create Targets / Assays**
 - **Edit Targets / Assays**
 - **Delete Targets / Assays**
- Analysis settings library
 - **Create Analysis Settings**
 - **Edit Analysis Settings**
 - **Delete Analysis Settings**
- Dye library
 - **Add New Custom Dye**
 - **Delete Custom Dye**

Apply a filter to search a library

You can filter the Sample, SNP assay, Target, and Analysis Settings Libraries.

1. Access a library of interest.
 - In the menu bar, select **Tools ▶ (Library of interest)**.
 - In the **Plate** tab of an open EDT or EDS file, click **Advanced Setup**, then select  **Action ▶ Import from Library**.
2. Select a feature from the first dropdown list. Each column of the table is an available feature.

3. Select a condition from the second dropdown list to define the feature. The conditions will vary by feature.
4. Enter a value or text by which to filter.
5. Click **Apply Filter**.

Define samples in the Samples table

In the **Plate** tab, click **Advanced Setup**, then perform one of the following actions in the **Samples** table.

Option	Action
Manually define a sample	a. Click  Add. b. Click a cell to edit the attributes for the sample. c. (Optional) Click  in the table header to add a Custom Attribute column. a. Click the Custom Attribute column header, then edit the header with a new sample attribute. b. Click a Custom Attribute cell in the table, then enter the attribute information.
Import samples from a TXT or XLS file	a. Select  Action ▶ Import from File. b. Navigate to and select a file, then click Open.
Import samples from the Sample Library	a. Select  Action ▶ Import from Library. b. (Optional) Apply a filter to search for a specific sample (see page 111). c. Select one or more samples, then click Add Selected. Note: Shift-click or Ctrl-click to select multiple samples.
Save a sample to the Sample Library	Select a sample row, then select  Action ▶ Save to Library .
Delete a sample from the table	Select a sample row, then click  .



Define targets in the Targets table

In the **Plate** tab, click **Advanced Setup**, then perform one of the following actions in the **Targets** table.

Option	Action
Manually define a target	- Click  Add. - Click a cell to edit the attributes for the target.
Import targets from a TXT or XML file	- Select  Action ▶ Import from Library. - Click Import or Import AIF. - Navigate to and select a file, then click Import. - Select one or more targets, then click Add Selected. Note: Shift-click or Ctrl-click to select multiple targets.
Import targets from the Target Library	- Select  Action ▶ Import from Library. (Optional) Apply a filter to search for a specific target (see page 111). - Select one or more targets, then click Add Selected. Note: Shift-click or Ctrl-click to select multiple targets.
Save a target to the Target Library	Select a target row, then select  Action ▶ Save to Library .
Delete a target from the table	Select a target row, then click  .

Define SNP assays in the SNP Assays table

In the **Plate** tab, click **Advanced Setup**, then perform one of the following actions in the **SNP Assays** table.

Option	Action
Manually define a SNP assay	- Click  Add. - Click a cell to edit the attributes for the SNP assay.
Import SNP assays from a TXT or XML file	- Select  Action ▶ Import from Library. - Click Import or Import AIF. - Navigate to and select a file, then click Import. - Select one or more SNP assays, then click Add Selected. Note: Shift-click or Ctrl-click to select multiple SNP assays.
Import SNP assays from the SNP Assay Library	- Select  Action ▶ Import from Library. (Optional) Apply a filter to search for a specific SNP assay (see page 111). - Select one or more SNP assays, then click Add Selected. Note: Shift-click or Ctrl-click to select multiple SNP assays.
Save a SNP assay to the SNP Assay Library	Select a SNP assay row, then  Action ▶ Save to Library .
Delete a SNP assay from the table	Select a SNP assay row, then click  .



Configure analysis settings

■ Guidelines for the analysis settings	115
■ View and configure the analysis settings	116
■ Ct settings overview	116
■ Flag settings overview	117
■ Advanced settings overview	117

This section describes the analysis settings that apply to all experiment types, unless otherwise noted.

Guidelines for the analysis settings

- We recommend analyzing the experiment with the default analysis settings.
- If the default analysis settings are not suitable for the experiment, modify the settings in the **Analysis Settings** dialog box, then reanalyze the experiment.
- Save modified analysis settings to the **Analysis Settings Library**.

The default analysis settings are different for each experiment type. The analysis settings determine for following parameters.

- How the baseline, threshold, and threshold cycle (C_t) are calculated
- Which flags are enabled
- Other analysis options that are specific to an experiment type

For detailed information about different types of analysis settings, see the following sections.

- “Ct settings overview” on page 116
- “Flag settings overview” on page 117
- “Advanced settings overview” on page 117
- “Standard curve settings overview” on page 52
- “Relative quantification settings overview” on page 65
- “Call settings overview (genotyping)” on page 71
- “Call settings overview (presence / absence)” on page 78
- “Melt curve settings overview” on page 83

View and configure the analysis settings

If security is enabled in the software, your account must have the permission of **Edit Analysis Settings** to configure the settings.

1. In the top-right corner of the **Results** tab, click .
2. View and (optional) configure the analysis settings.
For step-by-step instructions for adjusting analysis settings, see the desktop software Help.
3. Click **Apply**.
4. Click **Analyze** to reanalyze to experiment with the new settings.
5. (Optional) To save the settings in the **Analysis Settings Library**, click **Save**.
6. (Optional) To return to the default settings, click **Revert**.

C_t settings overview

The default C_t settings are appropriate for most applications. Configuration of the settings is an option for analysis of atypical or unexpected run data.

For step-by-step instructions for adjusting the C_t settings, see the desktop software Help.

Note: The **C_t Settings** feature is not available for experiments without a PCR stage, such as melt curve experiments.

Table 32 C_t Settings

Setting	Description
Data Step Selection	Determines the stage/step combination for C _t analysis (when there is more than one data collection point in the run method).
Algorithm Settings – Baseline Threshold	The Baseline Threshold Algorithm is used to calculate the C _t values. This algorithm is an expression estimation algorithm that subtracts a baseline component and sets a fluorescence threshold in the exponential region.
Algorithm Settings – Relative Threshold	The Relative Threshold (C _{rt}) Algorithm is used to calculate the C _{rt} values.
Default C_t Settings	Determines how the Baseline Threshold Algorithm is set. The Default C_t Settings are used for targets unless they have custom settings. For recommendations on adjusting baseline and threshold settings, see Table 33.

Table 32 Ct Settings (continued)

Setting	Description
C _t Settings for Target	• Default Settings selected—The Default C_t Settings are used to calculate the C_t values for the target. • Default Settings deselected—The software allows manual setting of the baseline or the threshold. For recommendations for adjusting baseline and threshold settings, see Table 33.

Table 33 Recommendations for manual threshold and baseline settings

Setting	Recommendation
Threshold	Enter a value for the threshold so that the threshold is: • Above the background. • Below the plateau and linear phases of the amplification curve. • Within the exponential phase of the amplification curve.
Baseline	Select the Start Cycle and End Cycle values so that the baseline ends before significant fluorescence signal is detected.

Flag settings overview

Use the **Flag Settings** to configure the following parameters.

- Adjust the sensitivity so that more wells or fewer wells are flagged.
- Change the flags that are applied by the software for each experiment type.

For step-by-step instructions for configuring the flag settings, see the desktop software Help.

Advanced settings overview

Use the **Advanced Settings** tab to change baseline settings for individual wells.

For step-by-step instructions for adjusting the advanced settings, see the desktop software Help.

Note: The **Advanced Settings** feature is not available for experiments without a PCR stage, such as melt curve experiments.



Documentation and support

Related documentation for the QuantStudio™ 1 Real-Time PCR System

Document	Publication number
<i>QuantStudio™ 1 Real-Time PCR System Installation, Use, and Maintenance Guide</i>	MAN0017853
<i>QuantStudio™ Design and Analysis Desktop Software Command-Line Application Guide</i>	MAN0010409
<i>QuantStudio™ 1 Real-Time PCR System Site Preparation Guide</i>	MAN0017854

Related documentation for the QuantStudio™ 3 Real-Time PCR System and QuantStudio™ 5 Real-Time PCR System

Document	Publication number
<i>QuantStudio™ 3 and 5 Real-Time PCR Systems Installation, Use, and Maintenance Guide</i>	MAN0010407
<i>QuantStudio™ 5 Real-Time PCR System Installation, Use, and Maintenance Guide (for China)</i>	MAN0028412
<i>QuantStudio™ Design and Analysis Desktop Software Command-Line Application Guide</i>	MAN0010409
<i>SAE Administrator Console v2.1 User Guide</i>	MAN0028414
<i>QuantStudio™ 3 and 5 Real-Time PCR Systems Site Preparation Guide</i>	MAN0010405

Customer and technical support

Visit thermofisher.com/support for the latest service and support information.

- Worldwide contact telephone numbers
- Product support information
 - Product FAQs
 - Software, patches, and updates
 - Training for many applications and instruments

- Order and web support
- Product documentation
 - User guides, manuals, and protocols
 - Certificates of Analysis
 - Safety Data Sheets (SDSs; also known as MSDSs)

Note: For SDSs for reagents and chemicals from other manufacturers, contact the manufacturer.

Limited product warranty

Life Technologies Corporation and/or its affiliate(s) warrant their products as set forth in the Life Technologies' General Terms and Conditions of Sale at www.thermofisher.com/us/en/home/global/terms-and-conditions.html. If you have any questions, please contact Life Technologies at www.thermofisher.com/support.

Glossary

biological replicates

Reactions that contain identical components and volumes, but evaluate separate samples of the same biological source (for example, samples from three different mice of the same strain, or separate extractions of the same cell line or tissue sample).

For runs that use biological replicate groups in a gene expression project, the values displayed in the biological replicates lists are calculated by combining the results of the separate biological samples and treating this collection as a single population (that is, as one sample).

For ΔC_t computations (normalizing by the endogenous control) in a singleplex experiment, the software averages technical replicates. The averages from the technical replicates are then averaged together to determine the value for that biological replicate.

endogenous control

A gene that is used to normalize template differences and sample-to-sample or run-to-run variation.

endpoint read

See post-PCR read.

post-PCR read

In genotyping and presence/absence experiments, the part of the instrument run that occurs after amplification. In genotyping experiments, fluorescence data collected during the post-PCR read are displayed in the allelic discrimination plot and used to make allele calls. In presence/absence experiments, fluorescence data collected during the post-PCR read are displayed in the presence/absence plot and used to make detection calls. Also called endpoint read.

pre-PCR read

In genotyping and presence/absence experiments, the part of the instrument run that occurs before amplification. The pre-PCR read is optional but recommended. Fluorescence data collected during the pre-PCR read can be used to normalize fluorescence data collected during the post-PCR read.

reference sample

In relative standard curve and comparative C_t ($\Delta\Delta C_t$) experiments, the sample used as the basis for relative quantification results. Also called the calibrator.

target

The nucleic acid sequence that is amplified and detected during PCR.

task

In the software, the type of reaction performed in the well for the target.

technical replicates

Reactions that contain identical components and volumes, and that evaluate the same sample; important for evaluating precision.

Index

A

- absolute quantity 46
- accessories 23
- Advanced Setup 107
- advanced temperature zones 104
- algorithm settings
 - Baseline Threshold Algorithm 116
 - Relative Threshold Algorithm 116
- Allelic Discrimination Plot, overview 69
- Amplification Plot
 - baseline 32
 - baseline setting 34
 - exponential phase 32
 - linear phase 32
 - optimize to view negative controls 36
 - overview 32
 - plateau phase 32
 - shape 32
 - threshold settings 33
- analysis settings
 - advanced settings 116, 117
 - configure 115
 - flag settings 116, 117
 - guidelines 115
 - melt curve settings 83
 - standard curve settings 116
- Analysis Settings Library 111
- Application 85
- assign
 - biological replicate group 110
 - sample 20, 107
 - SNP assay 20, 107
 - target 20, 107
 - task 20, 107, 108
 - well details 107
- assign standard dilutions 110
- assign targets, samples, and biological replicates,
from XLS file 101
- Audit 85
- Auto Delta, settings 104
- auto-analysis 31
- automatic threshold 116

B

- barcode scan 18
- Baseline Threshold Algorithm 116
- biological replicate groups, define and assign 110

C

- calibration data, override 42
- calibration results, view 41
- column in XLS file, well number 101
- comparative Ct experiments
 - compared to relative standard curve 56
 - endogenous control 55
 - Gene Expression Plot 63
 - negative controls 55
 - overview and purpose 55
 - PCR options 56
 - QC Plot 64
 - reaction types 55
 - reference genes 55
 - relative quantification settings 65
 - replicates 55
 - review results workflow 60
 - set up in software 58
- consumables 23
- Cq vs Ct vs Crt 28
- Ct precision, improve 35
- Ct settings
 - algorithm settings 116
 - default 116
 - select data step 116
 - tab 116
 - target 116
- custom dye 19
- custom experiments
 - overview 97
 - set up in software 98

D

- define
 - biological replicate group 110
 - samples 112
 - SNP assays 114

- targets 113
- delete files 93
- documentation, related 118
- Dye Library 111
- dye signal accuracy, Multicomponent Plot 38

E

- E-signature 85
- endogenous control
 - comparative Ct experiments 64
 - relative standard curve experiments 64
- example templates 102
- experiment
 - export 42, 45
 - setup workflow 16
- experiment data, import 100
- experiment template 102
- experiment types 11
- experiments, define parameters 21
- export
 - experiment or results 45
 - experiments or results 42

F

- flags 40

G

- Gene Expression Plot
 - comparative Ct experiments 63
 - overview 63
 - relative standard curve experiments 63
 - RQ vs Sample 63
 - RQ vs Target 63
- genotyping experiments
 - Allelic Discrimination Plot 69
 - call settings 71
 - manual call 71
 - PCR options 67
 - Presence/Absence Plot 76
 - reaction types 67
 - review results workflow 69
 - set up in software 67

I

- identify well problems 37
- Index Term 14
- instrument, power on 93

K

- kits 23

L

- library
 - apply filter 111
 - list of 107
- limited product warranty 119
- load plate 94, 95
- locked template, save 22

M

- melt curve experiments
 - overview and purpose 80
 - ramp increment 106
 - reaction types 80
 - reactions 82
 - review results workflow 82
 - set up in software 81
 - settings 83
- Melt Curve Plot
 - overview 82
 - view and assess 83
- method, create or edit 103
- monitor run 24
- Multicomponent Plot
 - dye signal accuracy 38
 - view 38

N

- network configuration and security 9
- NULL values in well assignments 101

O

- omit wells 35
- optical filter selection 19
- optical filter settings 105
- Overview 85

P

- parts and materials 23
- password security 9
- pause, add or adjust 105
- PCR options
 - comparative Ct experiments 56
 - genotyping experiments 67
 - presence/absence experiments 74
 - relative standard curve experiments 54

- standard curve experiments 47
- PCR reactions
 - guidelines for set up 24
 - handle samples and reagents 23
 - prepare 23
 - set up and run 48, 59, 68, 75
- permissions, role 88
- plate set up 20
- plate setup data, import 100
- Plate tab 107
- post-run summary 26
- PowerPoint 44
- preferences, set 15
- presence/absence experiments
 - call settings 78
 - IPC 73
 - overview and purpose 73
 - PCR options 74
 - reaction types 73
 - review results workflow 76
 - set up in software 74
- Presence/Absence Plot
 - overview 76
 - view and assess 77
- print run data and results 44

Q

- QC Plot
 - comparative Ct experiments 64
 - overview 64
 - relative standard curve experiments 64
 - view and assess 64
- QC Summary
 - review flags 40
 - view 40
- quantification cycle (Cq) 28
- Quick Setup 20, 107

R

- Raw Data Plot
 - determine signal accuracy 39
 - view 40
- reaction types
 - comparative Ct experiments 55
 - genotyping experiments 67
 - melt curve experiments 80
 - presence/absence experiments 73
 - relative standard curve experiments 54
 - standard curve experiments 46
- reagents, enter information 18
- real-time plots
 - view on instrument touchscreen 26

- zoom 26
- related documentation 118
- relative quantification settings 65
- relative quantity 53, 55
- relative standard curve experiments
 - assign standard dilutions 110
 - compared to comparative Ct 56
 - endogenous control 54
 - Gene Expression Plot 63
 - overview and purpose 53
 - PCR options 54
 - QC Plot 64
 - reaction types 54
 - relative quantification settings 65
 - replicates 54
 - review results workflow 60
 - set up in software 57
 - standard dilution series 54
- Relative Threshold Algorithm 116
- required components
 - negative controls 46, 54, 67
 - reference genes 54
 - replicates 46, 67
 - sample types 46, 67
 - standard 46
 - standard dilution series 46
- results, export 42, 45
- Results tab 29
- review
 - Amplification Plot 32
 - Multicomponent Plot 38
 - QC Summary 40
 - Raw Data Plot 39
 - Standard Curve Plot 50, 61
 - Well Table 37
- role
 - default 88
 - permissions 88
- run data, print 44
- run history 93
- run results
 - print 44
 - review workflow 31

S

- sample
 - assign to well 20, 107
 - define 112
- sample definition file
 - create 99
 - spreadsheet 99
- sample information, import 99
- Sample Library 111

- Samples table, manage 112
- save template 22
- Security 85
- Security, Auditing, E-signature Administrator Console, Overview 85
- select wells 21
- signal accuracy, Raw Data Plot 39
- SNP assay
 - assign to well 20, 107
 - define 114
- SNP Assay Library 111
- SNP Assays table, manage 114
- standard 46, 54
- standard curve, use a different standard curve 52, 63
- standard curve experiments
 - assign standard dilutions 110
 - overview and purpose 46
 - PCR options 47
 - reaction types 46
 - review results workflow 49
 - set up in software 47
- Standard Curve Plot
 - overview 50, 61
 - view and assess 51, 62
- standard curve settings 52, 63
- standard dilutions, set up 109
- standards, define and set up 109
- start run 24
- support, customer and technical 118

T

- target
 - assign to well 20, 107
 - define 113
- Target Library 111
- Targets table, manage 113
- task
 - assign to well 20, 107, 108

- descriptions 108
- template
 - adjust reaction volume 19
 - create 17
 - enter properties 17
 - information saved 102
 - open 17
 - optical filter selection 19
 - set up 16
 - use default method 19
- terms and conditions 119
- thermal protocol, create or edit 103
- transfer files 93

U

- unload plate 94, 95
- unlocked template, save 22

V

- VeriFlex, settings 104

W

- warranty 119
- well details
 - assign 107
 - define 107
 - view on instrument touchscreen 26
- Well Table
 - group or sort 37
 - overview 37
- wells, select 21
- workflow
 - review run results 31
 - run experiment 16
 - set up experiment 16

Z

- zoom 26

