



Anaesthesia as advanced as your MRI

FABIUS® MRI

The MRI challenge



Flexibility is a great thing. In the proper setting, flexibility can help save both time and money, helping to improve overall efficiency. Every now and then, however, the need for a dedicated solution arises. Some applications are so demanding, so specialized, that they require a solution that was designed solely to serve that special purpose.

MRI is one of those specialized tasks. Restricted access, high throughput, limited space and strong magnetic fields can pose challenging workflow problems. In many parts of the world, hybrid diagnostic/therapeutic suites are emerging as the new gold standard for advanced treatment strategies. When one considers the high cost and complex logistics of building a modern MRI suite, it quickly becomes obvious that there is little or no margin for error. Once you build it, it's pretty much set in stone. It has to be done right the first time.

ANAESTHESIA AND MRI

There are many reasons why anaesthesia plays such a pivotal role in MRI diagnostics. A significant number of patients are critically ill and require ventilation support during both transport and imaging. Small children and neonates generally require sedation in order to eliminate motion artefacts. Neurosurgeons are now using MRI directly in the operating room to provide integrated planning and guidance for new and complex procedures. Other patients, while otherwise stable, suffer from claustrophobia and can't tolerate the confined conditions inside a tomography unit without sedation. There's no doubt: to get the most out of any MRI diagnostic unit, you need integrated anaesthesia. But the question remains: Is your anaesthesia workstation up to the challenge of the MRI environment?

ADVANCED ANAESTHESIA


FOR ADVANCED APPLICATIONS: THE FABIUS MRI

The Fabius MRI features our high-precision E-Vent piston ventilator technology, renowned for both performance and safety. Its electronically controlled, electrically driven ventilator requires no drive gas, making it particularly well suited for the typical MRI environment. The Fabius MRI puts all the major ventilation modes at your fingertips, from volume and pressure controlled, pressure support and SIMV/PS. It lets you ventilate with ICU-like performance, providing enhanced safety, confidence and control and is equally suited for adult, paediatric and neonatal patients. With advanced features such as dynamic compliance compensation, fresh gas decoupling and a pop-up APL release valve, you can rest assured that ventilation is gentle, safe and precise. With its leak-tight breathing system, you can even use the Fabius MRI for low-flow anaesthesia. But that's just the beginning. To make the Fabius MRI the ventilator of choice for MRI applications, it's been engineered and built especially for use in powerful magnetic fields. In the fact, the MR conditional Fabius MRI has been tested for use with MRI systems of 1.5 and 3 Tesla and can be safely operated at field strengths of 40 mtesla (400 gauss) from the tomography unit. This puts an end to extra-long breathing circuits with all the compliance and dead space issues involved. Now you can perform valuable diagnostics with enhanced confidence and safety.



MADE WITH ERGONOMICS IN MIND

The Fabius MRI doesn't just provide excellent ventilation therapy. It was also designed to enhance your workflow. A large, high-visibility colour screen keeps you informed on patient and ventilator status even from a distance. Integrated warning and alarm LEDs provide additional visual support. The intuitive menu structure and the familiar, 3-step Dräger operating philosophy – select, adjust, confirm – will move you almost effortlessly through its systems. Three large drawers built in to the trolley provide ample storage space for additional equipment and consumables. The external COSY breathing system features a simple design with very few parts for easy assembly and cleaning. It can be mounted on either side to precisely complement your individual needs. Equipped with our standard CLIC absorber, you can change soda lime canisters quickly and cleanly with no interruption of operation.

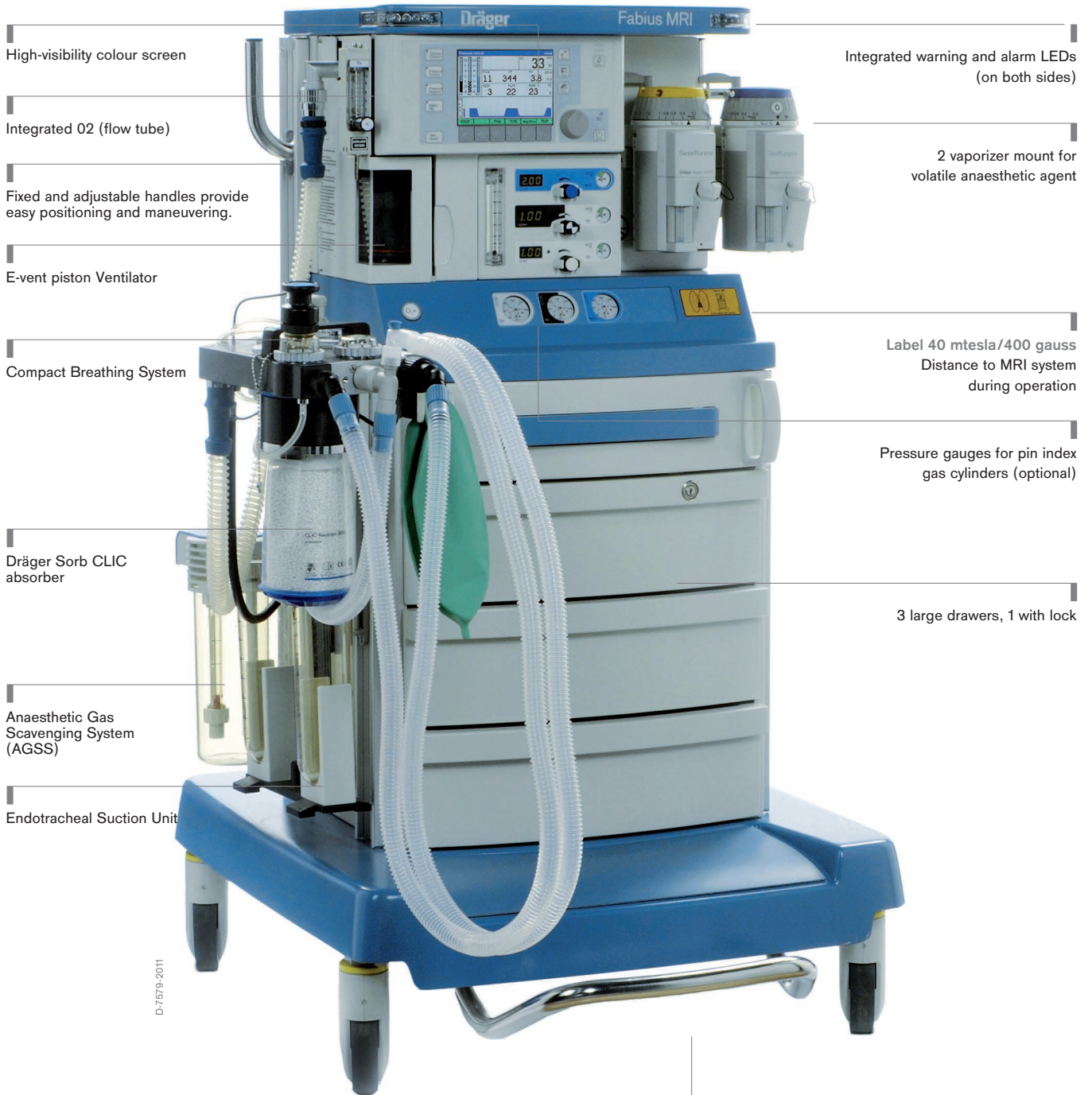

High-visibility colour screen

Integrated O2 (flow tube)

Fixed and adjustable handles provide easy positioning and maneuvering.

E-vent piston Ventilator

Compact Breathing System

Dräger Sorb CLIC absorber

Anaesthetic Gas Scavenging System (AGSS)

Endotracheal Suction Unit

Integrated warning and alarm LEDs (on both sides)

2 vaporizer mount for volatile anaesthetic agent

Label 40 mtesla/400 gauss
Distance to MRI system during operation

Pressure gauges for pin index gas cylinders (optional)

3 large drawers, 1 with lock

D-7579-2011

Central brake effectively enables fast and easy interlock or release for enhanced safety especially in magnetic fields.



COMPACT BREATHING SYSTEM „COSY 2.6“

- Clic-Adapter ready (Single use absorber)
- Short, flexible COSY mounting arm for enhanced ergonomics
- Integrated cable management
- Left- and right side mounting possible for optimal adaption and support of your MRI environment.



HIGH RESOLUTION TFT COLOR DISPLAY (6.5" DIAGONAL)

- Enhanced visibility thanks to optimized brightness and contrast
- A single screen controls all functions and monitors ventilation
- Oxygen monitoring
- Toggle between control and curve screen
- Standard Dräger user interface and operation



2 POWERFUL ADDITIONAL LEDS INTEGRATED IN THE TOP PLATE

- WARNING (yellow-flashing) and
- ALARM (red-flashing)
- Visibility from different angles and from a distance
e.g. the control room



MRI Suite and control room

The cockpit's view from the operator control area. Everything in view and under control.

WE UNDERSTAND

Over the years, Dräger products have made such an impact on hospital design that it's hard to imagine life without them. The name Dräger stands for quality and innovation the world over. So when you decide it's time to add an MRI suite at your hospital, consider Dräger. We understand the special challenges involved in designing and building an MRI environment. Critical factors such as sufficient shielding, adequate access, ferro-magnetic materials and the need for specialized equipment – anaesthesia workstations, monitors and even light bulbs – all have to be taken into account. Architectural considerations for the control room, patient holding area, changing rooms and of course the tomography unit itself have to be carefully thought through in order to achieve the best possible workflow for any given situation. The goal is to create an environment which lets you work smoothly and efficiently in spite of all the restrictions, helping you to get the most out of your suite.

NEW FRONTIERS

Modern advances in medical technology have brought about changes in medical philosophy. In many disciplines, the line between diagnostics and therapy is slowly disappearing. So-called "hybrid suites" or "hybrid theatres" are beginning to emerge as the method of choice to provide integrated high-resolution imaging and precision therapy in a single package. Here too, Dräger can assist you with both the technology and the know-how to create just such a solution for your hospital.

Anaesthesia as advanced




Induction Room
Induction of anaesthesia in a patient before entering the Neurosurgical MRI-OR

Patient transport
Transport patient into the MRI-OR

SPACE PLANNING IN THE MRI ENVIRONMENT

Today, the framework for anaesthesia in the MRI environment is quite different from that found in the operating room. Anaesthesia management for MRI doesn't just involve simply duplicating comparable operating room requirements. Instead, its workflow design includes a dedicated space for patient pathways (this applies in particular to intensive care patients) which allows for proper patient care depending on the individual conditions and challenges involved. For this reason, a comparable space reserve and at least the device assembly space available in an operating room area have to be taken into account. In addition to the separate areas for patient preparation and recovery, the design of the transport path to the MRI must be barrier-free. The most direct access possible, as well as the paths between the rooms (e.g. preparation room, MRT room, recovery), which are generally separated by sliding doors, must be designed to accommodate patients on transport stretchers.

Further information regarding space planning and planning support can be found in the separate product information brochure for space planning in MRI environments.

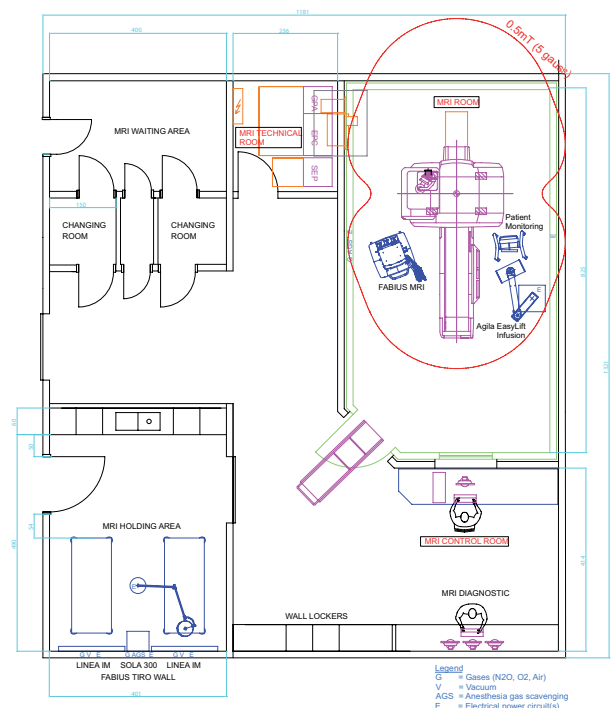

as your MRI



Neurosurgical MRI-OR
 Anaesthesia maintenance phase. Patient monitoring during MR sequences
 (Anaesthesia Device, Patient Monitor, Syringe Pumps)



Anaesthesia maintenance phase
 Fabius MRI in use during neurosurgery procedure.

DRAWING ON A WORLD OF KNOWLEDGE

To do this, we not only rely on our own extensive experience. We also draw on the expertise of strong partners who specialize in patient monitoring and infusion technology for MRI environments. We work together closely with architects and contractors to develop a highly efficient and ergonomic design that is best suited to your individual needs.

By taking advantage of a dedicated MRI solution, you can increase your overall efficiency by reducing transfer and setup times in the MR room. The time you save here can be put to better use – performing diagnostics.

PUTTING IT ALL TOGETHER

Let Dräger help you plan, build and run the MRI solution that's right for you. It's an investment in the future you won't regret. Our representatives will be happy to show you the most efficient and effective answer to your MRI challenge.





**NOW MAKE THE RIGHT CHOICE FOR SERVICE:
DRÄGERSERVICE®**

DrägerService offers a wide range of flexible service options which you can select and combine to meet your individual needs – from technical support to complete, fixed-budget service contracts. Remote Service solutions offer additional potential to improve the uptime of devices.

Accessories & Consumables

A wide array of MRI conditional accessories and consumables are available to choose from, giving you the possibility to tailor the Fabius® MRI to create exactly the anaesthesia solution you need:



D-559-2009

VentStar® MRI 300

Disposable breathing circuit, consisting of 2 smoothbore hoses, Y-piece and LuerLock elbow.

Suitable for use in MRI environment.

Length: 300 cm (118 inch).

Latex free.



D-34634-2009

VentStar® MRI (N) 300

Disposable neonatal breathing circuit, consisting of 2 smoothbore hoses (Ø 10 mm), angled Y-piece with LuerLock.

Suitable for use in MRI environment.

Length: 300 cm (118 inch).

Latex free.



MT-3716-2008

VacuSmart® Gel

The VacuSmart Gel is a disposable cartridge with an integrated gelling agent for all Dräger Medical bronchial aspirators using 700 ml secretion jars. The gelling agent coagulates the bronchial secretion and prevents thereby the leakage e.g. in a case of disposal in a waste press.



MT-2667-2004

Spirolog® Flowsensor

Spirolog Flowsensor is a hot-wire sensor for measuring volumetric gas flow delivered by Dräger anaesthesia devices and ventilators.



MT-2004-2008

CLIC Absorber 800+

Disposable CLIC Absorber 800+ is filled with Drägersorb 800+, 1.2 l (42.2 fl oz.).

CORPORATE HEADQUARTERS

Drägerwerk AG & Co. KGaA
Moislinger Allee 53–55
23558 Lübeck, Germany

www.draeger.com

Manufacturer:

Dräger Medical GmbH
Moislinger Allee 53–55
23558 Lübeck, Germany

As of August 2015:

Dräger Medical GmbH changes
to Drägerwerk AG & Co. KGaA.

REGION EUROPE CENTRAL AND EUROPE NORTH

Dräger Medical GmbH
Moislinger Allee 53–55
23558 Lübeck, Germany
Tel +49 451 882 0
Fax +49 451 882 2080
info@draeger.com

REGION EUROPE SOUTH

Dräger Médical S.A.S.
Parc de Haute Technologie d'Antony 2
25, rue Georges Besse
92182 Antony Cedex, France
Tel +33 1 46 11 56 00
Fax +33 1 40 96 97 20
dlmfr-contact@draeger.com

REGION MIDDLE EAST, AFRICA

Dräger Medical GmbH
Branch Office
P.O. Box 505108
Dubai, United Arab Emirates
Tel +971 4 4294 600
Fax +971 4 4294 699
contactuae@draeger.com

REGION ASIA / PACIFIC

Draeger Medical
South East Asia Pte Ltd.
25 International Business Park
#04-27/29 German Centre
Singapore 609916, Singapore
Tel +65 6572 4388
Fax +65 6572 4399
asia.pacific@draeger.com

REGION CENTRAL AND SOUTH AMERICA

Dräger Panama Comercial S. de R.L.
Complejo Business Park,
V tower, 10th floor
Panama City
Tel +507 377 9100
Fax +507 377 9130
contactcsa@draeger.com

CANADA

Draeger Medical Canada Inc.
2425 Skymark Avenue, Unit 1
Mississauga, Ontario, L4W 4Y6
Tel +1 905 212 6600
Toll-free +1 866 343 2273
Fax +1 905 212 6601
Canada.support@draeger.com