

EasyVAC[®] PLUS VACUUM REGULATORS

The **EasyVAC[®] PLUS** continuous suction vacuum regulators range is suitable in all medical suction applications in hospitals and nursing homes. The device is made of a strong techno-polymer body, with a quick I/O switch button, a suction adjustment Soft Grip inserts knob for an easy handling with a "Push & Lock" position system and a vacuum control gauge with three possible end-of-scale choices: -250 mbar, -600 mbar and -1000 mbar. A protective silicone cover and a techno-polymer support prevent gauge damages caused by possible shocks received during transport or use. The **EasyVAC[®] PLUS** vacuum regulators are manufactured with outlet threaded connection for screwing in the usual collection containers for suctioned fluids, or for direct connection to the **EasySAFE[®] PLUS** safety jar, through a specific quick release integrated connection. The use of the **EasySAFE[®] PLUS** safety jar is highly recommended to get a complete suction device able to guarantee total protection of the vacuum regulator and plant.

Standard Inlet connection

EasyVAC[®] PLUS 250 - 600 - 1000 ▶ ISO G. 1/4" M.

Maximum suction flow

EasyVAC[®] PLUS 250 ▶ 50 L/min. ±5 L/min. at -220 mbar

EasyVAC[®] PLUS 600 ▶ 72 L/min. ±5 L/min. at -550 mbar

EasyVAC[®] PLUS 1000 ▶ 115 L/min. ±5 L/min. at -950 mbar

Maximum available negative pressure effective value

EasyVAC[®] PLUS 250 ▶ -220 mbar ±20 mbar

EasyVAC[®] PLUS 600 ▶ -550 mbar ±20 mbar

EasyVAC[®] PLUS 1000 ▶ -950 mbar

Over depression protecting valve

EasyVAC[®] PLUS 250 ▶ Preset, with anti clogging system

EasyVAC[®] PLUS 600 ▶ Preset, with anti clogging system

EasyVAC[®] PLUS 1000 ▶ N/A

Vacuum gauge

EasyVAC[®] PLUS 250 ▶ 0 ÷ -250 mbar

EasyVAC[®] PLUS 600 ▶ 0 ÷ -600 mbar

EasyVAC[®] PLUS 1000 ▶ 0 ÷ -1000 mbar

Maximum depression supply value


EasyVAC[®] PLUS 250 - 600 - 1000 ▶ -950 mbar

Suction value adjustment system

EasyVAC[®] PLUS 250 - 600 - 1000 ▶ With membrane

I/O SWITCH

Quick push switch button.

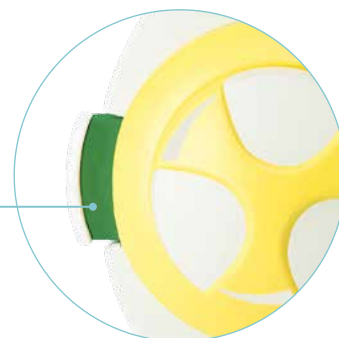

OUTLET CONNECTION

ISO G. 1/2" M. + quick release system for **EasySAFE[®] PLUS** safety jar connection only.

Vacuum gauge

EasyVAC[®] PLUS 250 ▶ 0.470 Kg

EasyVAC[®] PLUS 600 ▶ 0.370 Kg

EasyVAC[®] PLUS 1000 ▶ 0.370 Kg

Overall dimensions

safety jar excluded

EasyVAC[®] PLUS ▶ L 91 x W 106 x H 185 mm



VACUUM GAUGE READING SCALE

Vacuum gauge with colored sectors scale, for easy reading of set vacuum values.

VACUUM GAUGE PROTECTION

Silicon cover and anti shock support.

ADJUSTMENT ASSEMBLY

Quick release cover 1/12 turn.

VACUUM SETTING KNOB

With Push&Lock system and Soft Grip inserts for easy handling.



SAFETY JAR

EasySAFE® PLUS safety jar, optional.

The small sizes of the **EasyVAC® PLUS** vacuum regulators help in installing multiple sets of medical gas devices into the wall terminal units. The extreme rationality and simplicity, combined with the sophisticated technical execution, allow both operators and patients to appreciate the safety and reliability of this medical device.