


Item Code: 847

Indications

- Multifragmentary (dislocated) fractures proximal humerus
- Pseudarthrosis proximal humerus
- Osteotomies proximal humerus
- Non-union & mal-union proximal humerus

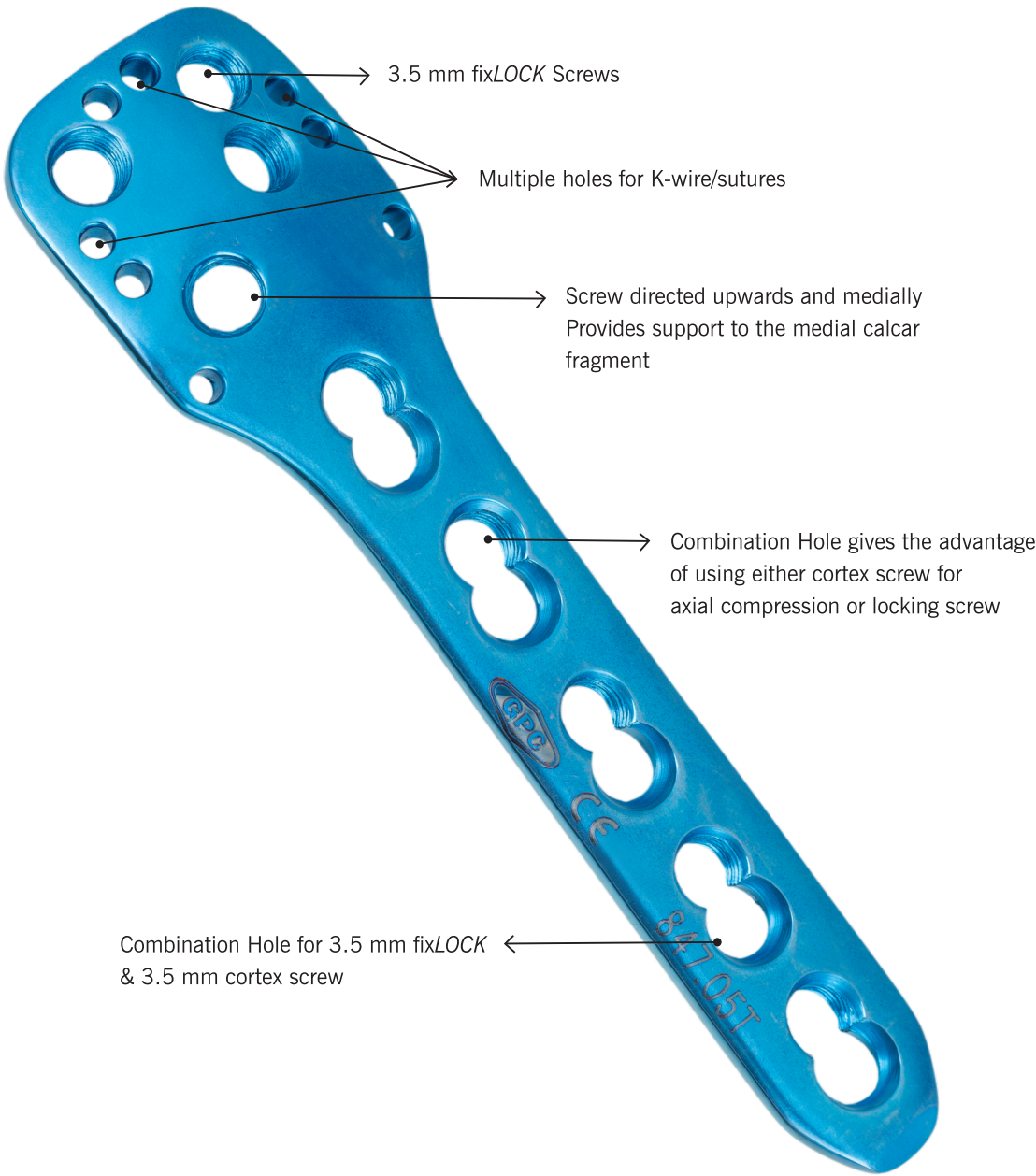

Features

- Universal for right & left humerus
- Low profile construct for minimum implant prominence
- Uniform stiffness of all segments, prolonging fatigue life of the implant
- Smaller area of post-fixation avascularity
- The combination of conventional and fixLOCK screws offers optimum fixation regardless of bone density
- Temporary reduction can be achieved by K-wire and the holes can later be utilized as suture holes
- Suture holes- For fixation of sutures or cerclage wires attached to rotator cuff muscles or the tubercles, it neutralizes the muscle tension and helps to maintain reduction
- Fixed angle stable construct with long oblique screw to stabilize medial calcar fragment
- Screw trajectory designed for optimal fixation of comminuted fractures
- Available in both Titanium & Stainless Steel

Sizes Available

No. of Holes	Code	
	S.S.	Titanium
3	847.03	847.03T
4	847.04	847.04T
5	847.05	847.05T
6	847.06	847.06T
8	847.08	847.08T
10	847.10	847.10T
12	847.12	847.12T