


**Elastica-Van Gieson kit**

Four-reagent kit for staining elastic fibers and differentiation between elastic tissue, collagen and other types of connective tissue. The rapid method enables a satisfactory result with shorter section staining time.


for 100 tests      WVGB-100T  
4x100 mL      WVGB-K-100

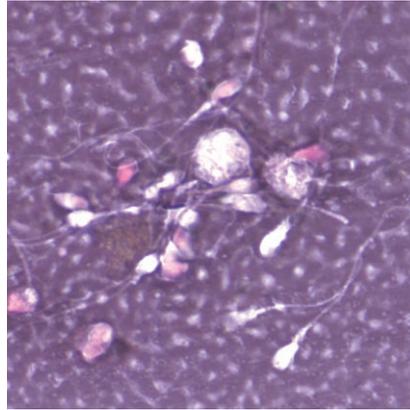

**Eosin and nigrosine vital kit**

Two-reagent kit containing separate dyes for rapid detection of sperm vitality and simple visualization of dead and living sperm cells. The nigrosin stain provides dark background for easier recognition of both viable and non-viable spermatozoa.


2x30 mL      ENVK-K-30  
2x100 mL      ENVK-K-100

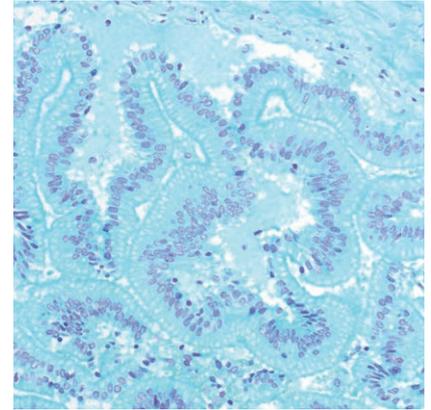

**Feulgen kit**

Five-reagent DNA staining kit according to Feulgen. For use in semiquantitative DNA determination in histological and cytological samples. Specimen is first treated with hydrochloric acid creating an aldehyde group of DNA that can be visualized by Schiff (BioSchiff) reagent. This reaction is specific for nuclear DNA.

for 100 tests      FE-100T  
6x100 mL      FE-K-100

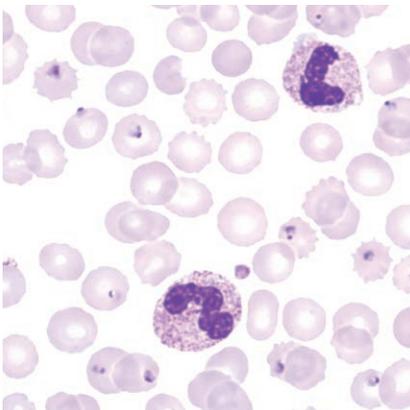

**Field kit**

Ready-to-use two-reagent kit for rapid and efficient staining and detection of parasites in hematology samples. Primarily used for staining thin and thick blood smears (dense drop) for purpose of diagnosing blood parasites. Reagents are stored in containers that can be used as staining jars.


for 100 tests      FD-100T  
2x500 mL      FD-K-500



**Fouchet-Van Gieson kit**

Three-reagent kit for visualization of bilirubin and collagen according to Kutlick. Bilirubin is a yellow-brown pigment, but changes to green due to oxidation induced by Fouchet reaction. Green bilirubin can easily be detected on yellowish and pink colored background.

for 100 tests      FVG-100T  
3x100 mL      FVG-K-100

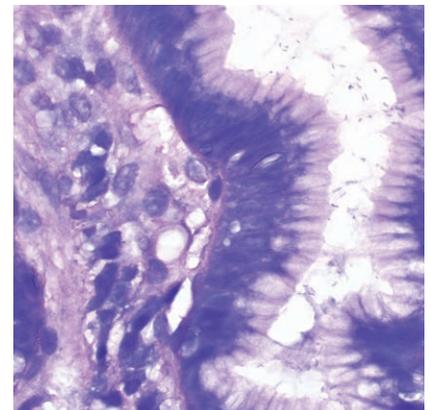

**Giemsa HP kit**

Four-reagent kit for staining *Helicobacter pylori* in gastroscopic sections according to Lennart. Advantages of this method for detecting *H. pylori* are sensitive and reproducible results and easy performance.

for 100 tests      GMHP-100T

