


OBTURA - TECHNICAL SPECIFICATIONS

Parameters	Specifications
Device Size (Fr)	6Fr and 8Fr
Device Total Length (mm)	205 ±10
Device Effective Length (mm)	155 ±10
Complete Degradation Period	90 Days
Anchor Material	Poly(D,L-lactide-co-glycolide)
Device Material	PDLG anchor, Poly Glycolic Acid (PGA) Suture and Bovine type-I Collagen plug
Device Size Compatibility	6Fr device compatible with 6Fr & 7Fr Introducer Sheath 8Fr device compatible with 8F & 9Fr Introducer Sheath
Introducer Sheath Compatibility	Medtronic: InPut® Introducer Sheath Cordis Corporation: AVANTI™ Introducer Sheath Terumo Corporation: PINNACLE® Introducer Sheath Oscor Inc: Adelante® Merit Medical: Prelude®

OBTURA - ORDERING INFORMATION

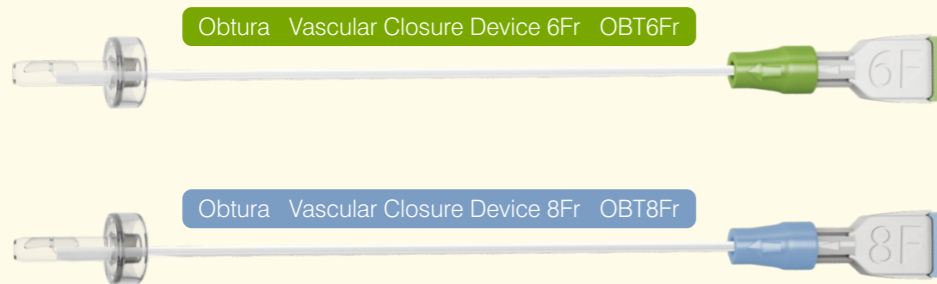

More to Life

Manufacturer:

Meril Life Sciences Pvt. Ltd.
Survey No. 135/139, Bilakhia House,
Muktanand Marg, Chala, Vapi - 396 191.
Gujarat, India.
T +91 260 240 8000

Subsidiary companies:

Meril Life Sciences Pvt. Ltd.
301, A-Wing, Business Square,
Chakala, Andheri Kurla Road,
Andheri East, Mumbai 400 093
T +91 22 39350700
F +91 22 39350777

Meril, Inc.

2436 Emrick Boulevard,
Bethlehem, PA 18020
T +610 500 2080
F +610 317 1672

Meril South America

Meril South America
Doc Med LTDA
Alameda dos Maracatins 992,
Edifício Aquarius, Andar 8,
Bloco B, Moema - CEP -04089001,
Sao Paulo, Brazil.
Tel: +5511 2506-3624

Meril GmbH.

Bornheimer Strasse 135-137,
D-53119 Bonn.
Germany.
T +49 228 7100 4000
F +49 228 7100 4001

Meril Tibbi Cihazlar

Meril Tibbi Cihazlar İmalat ve Ticaret A.Ş.
İçerenköy Mah. Çetinkaya Sok.
Prestij Plaza No:28
Kat:4 Ataşehir, 34752
İstanbul / Turkey
T +90 216 641 44 24
F +90 216 641 44 25

Meril China Co. Ltd.

2301b,23f, Lixin Plaza,
no 90,South Hubin Road,
Xiamen, China
T 0086-592-5368505
F 0086-592-5368519

Meril SA Pty. Ltd.

102, 104, S101 and S102,
Boulevard West Office Park,
142 Western Service Road,
Erf 813 Woodmead Extension 17
Sandton, Johannesburg – 2191
South Africa
T +27 11 465-2049
F +27 86 471 7941

Meril Medical LLC.

Nauchnyi Proezd 19,
Moscow, Russia – 117 246.
Office - +7 495 772 7643

EU representative.

Obelis S.A.
Bd, General Wahis 53,
1030, Brussels, Belgium.
T +32 2 732 5954
F +32 2 732 6003
E mail@obelis.net


E askinfo@merillife.com
W www.merillife.com



More to Life

Obtura™

Vascular Closure Device

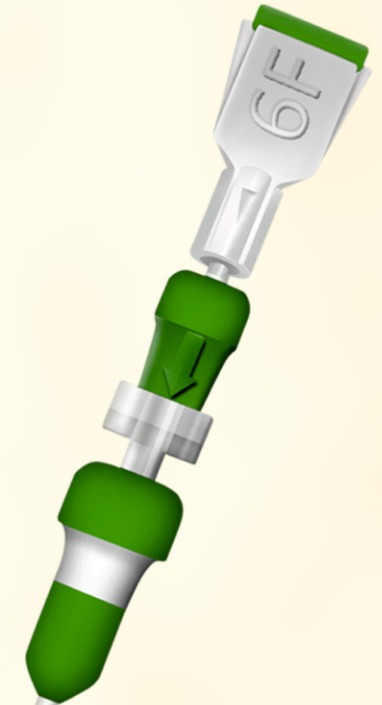

Contact your country local Meril sales representative for availability of sizes.
Obtura™ is a registered trademark of Meril Life Sciences Pvt. Ltd.
Obtura™ is not approved by USFDA and not available for sale in USA.



OBT/BROCHURE008/20211018/GLOBAL



Obtura™ is a collagen based convenient 3-step mechanical **Vascular Closure Device** for effective hemostasis.

No exchange of introducer sheath


Specially designed Obtura is compatible with the standard existing introducer sheath used during the interventional procedures leading to less blood loss

Rapid Hemostasis

Efficient Anchor-Collagen plug (Bovine type-I) sandwiches the arteriotomy and gives rapid hemostasis at puncture site


Rapid absorption

Obtura facilitates patient friendly sealing with the implant deployed at the puncture site which gets completely absorbed

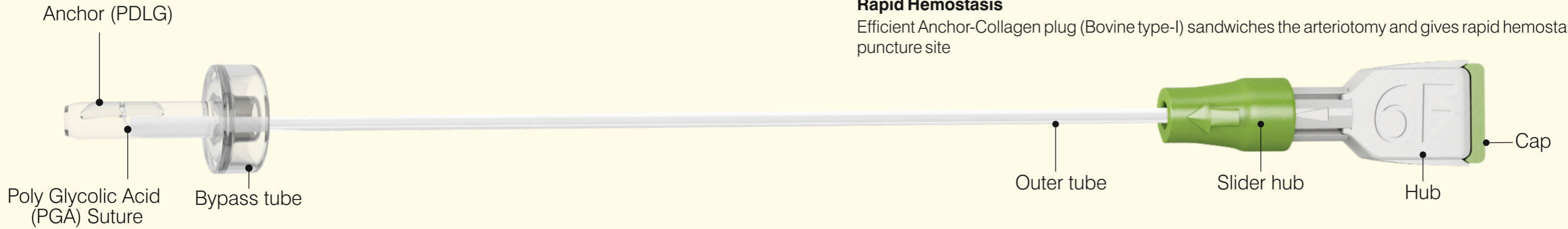


Unique 3-step deployment technique to reduce procedural time and blood loss


Deployment Result

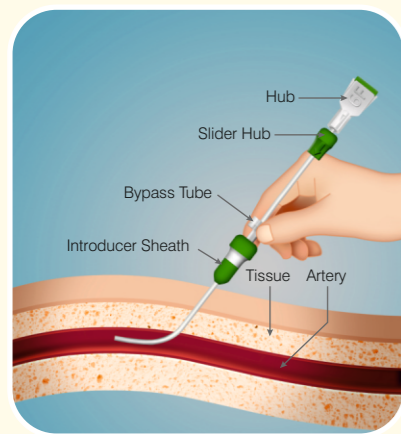

Step 1

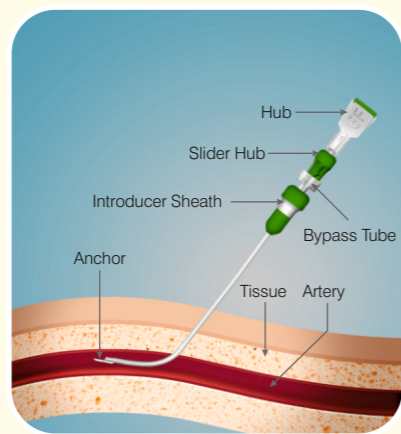

Step 2

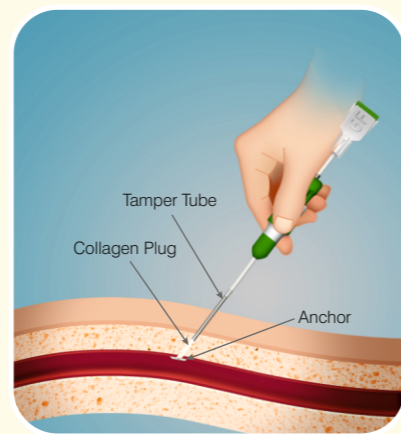

Step 3

