

A.80

BIL /BILE

Code: 0BIL
Code: 0BILE


Disposable needle for iliac-sternal aspiration
with anatomic grip (BIL) and
Explant version (BILE)

A



Anatomic grip

B

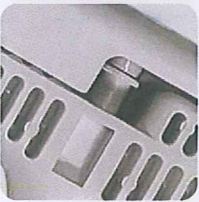


Luer Lock
connection

C



The threaded ring
nut can be inserted
or removed from
the handle

D

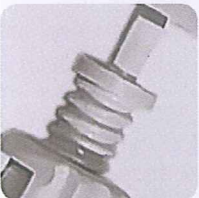


Millimetric feed

E




Triple sharpened tip

F

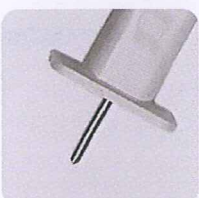


Depth limiter

G

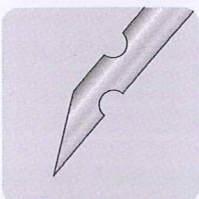

BILE version
tip with lateral
holes for explant

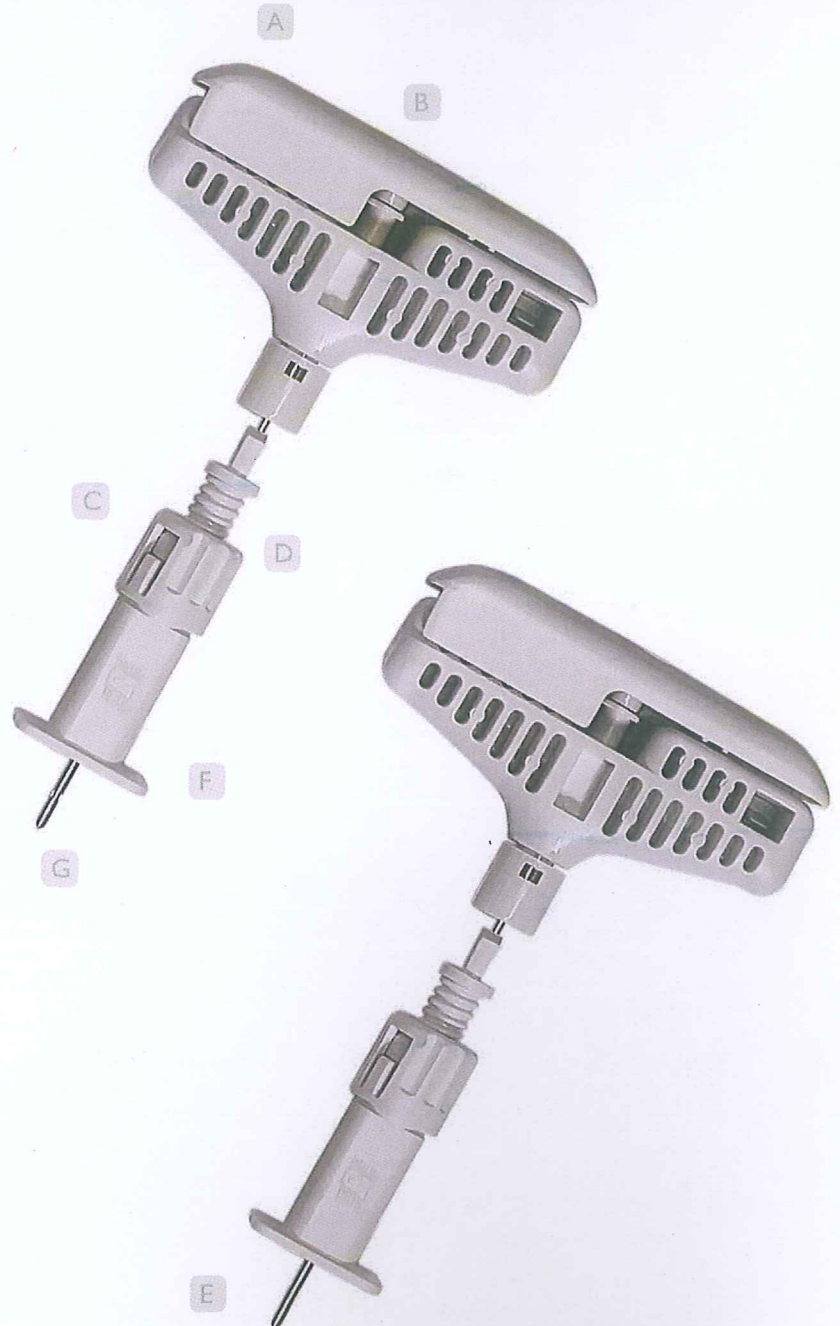


CE
0426

biomedical


BIL /BILE

Code: OBIL
Code: OBILE


Disposable needle for iliac-sternal aspiration
with anatomic grip (BIL) and
Explant version (BILE)

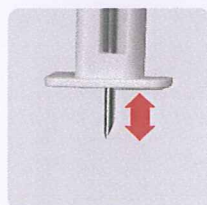


mm 17/32

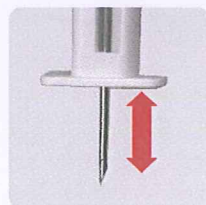

mm 36/51

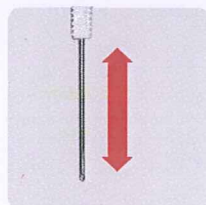

mm 55/70

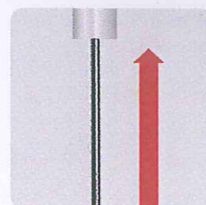

mm 80/95



BILE version


To order:

OBIL

CODE	SIZE	PACKAGE	C. CODE
OBIL 1417/080	14Gx17/80 mm	20	LIGHT GREEN
OBIL 1432/095	14Gx32/95 mm	20	LIGHT GREEN
OBIL 1517/080	15Gx17/80 mm	20	BLUE-GREY
OBIL 1532/095	15Gx32/95 mm	20	BLUE-GREY
OBIL 1617/080	16Gx17/80 mm	20	WHITE
OBIL 1632/095	16Gx32/95 mm	20	WHITE
OBIL 1817/080	18Gx17/80 mm	20	PINK
OBIL 1832/095	18Gx32/95 mm	20	PINK

OBILE

CODE	SIZE	PACKAGE	C. CODE
OBILE 1417/080	14Gx17/80 mm	20	LIGHT GREEN
OBILE 1432/095	14Gx32/95 mm	20	LIGHT GREEN
OBILE 1517/080	15Gx17/80 mm	20	BLUE-GREY
OBILE 1532/095	15Gx32/95 mm	20	BLUE-GREY
OBILE 1617/080	16Gx17/80 mm	20	WHITE
OBILE 1632/090	16Gx32/95 mm	20	WHITE
OBILE 1817/080	18Gx17/80 mm	20	PINK
OBILE 1832/095	18Gx32/95 mm	20	PINK

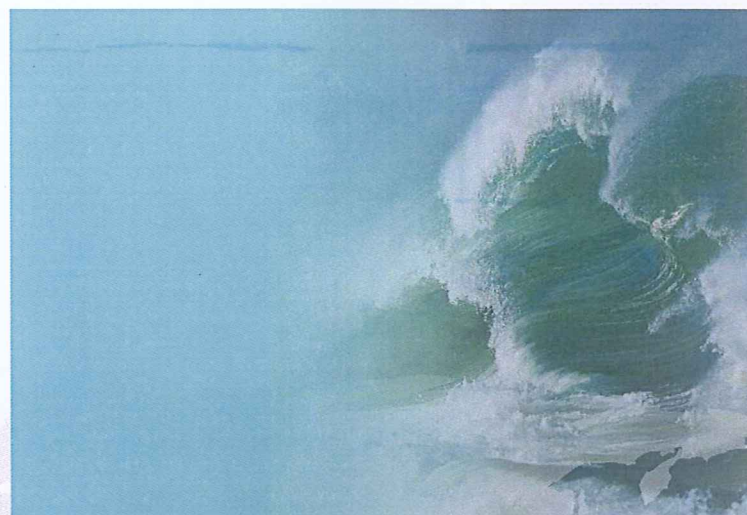

On request Biomedical is available to provide needles with different lengths.

biomedical

Biomedical srl Via G.B. Lulli, 43 50144 Firenze Italy
Tel +39 055 352678 +39 055 332822 Fax +39 055 331629
www.biomedical-bio.net sales@biomedical-bio.net

Bone biopsy needles

Guillotine needles for biopsy on soft tissues

Cyto-aspiration needles

Different specialized needles and sets


Loco-regional anaesthesia needles and sets

Thermic and absorbing sheets

Biocut-M

Reusable automatic device
for soft tissues biopsy

Code: Biocut-M

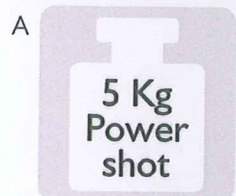

Charge
indicator
for cannula




Charge
indicator
for stylet



Shooting
button



Piston to
charge the
cannula

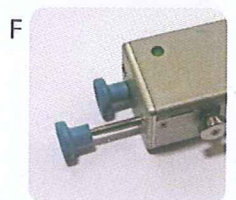

Piston to
charge the
stylet



Easy and safe

To be used with
PG needle



CE

biomedical



Biocut-M

Reusable automatic device
for soft tissues biopsy

Code: Biocut-M

Instructions of use:



Open the cover
and insert PG needle
into the gun



When the indicators
are green
Bioshot-M is charged




Close the cover

Withdraw the piston
on the left to charge
the cannula


PUSH



Push the shooting button

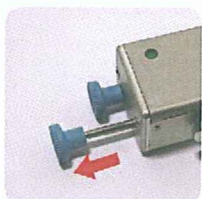

Withdraw the piston
on the right to
charge the stylet

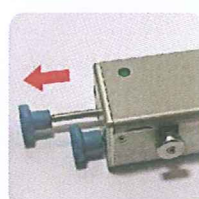


Withdraw
the piston on the left
to discharge
the sample

To order:

CODE
BIOCUT M

PACKAGE
I

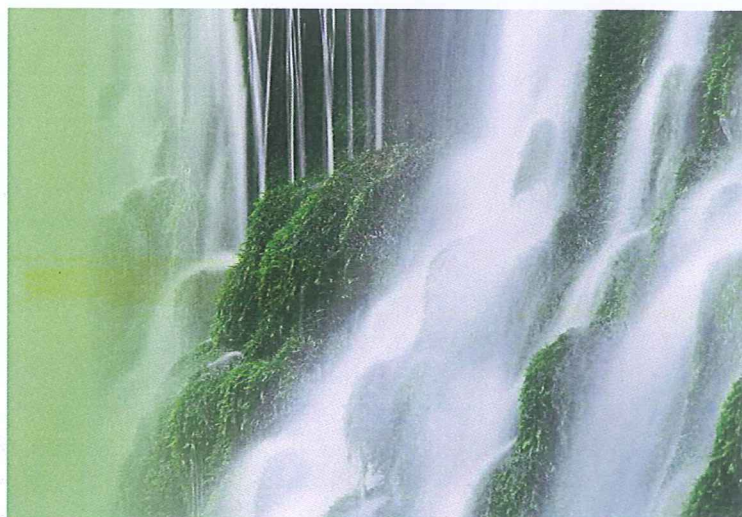

biomedical

Biomedical srl Via G.B. Lulli, 43 50144 Firenze Italy
Tel +39 055 352678 +39 055 332822 Fax +39 055 331629
www.biomedical-bio.net sales@biomedical-bio.net

Bone biopsy needles
Guillotine needles for biopsy on soft tissues
Cyto-aspiration needles
Different specialistic needles and sets
Loco-regional anaesthesia needles and sets
Thermic and absorbing sheets