



touchTymp

Our ALL-TOUCH Impedance Line



Simply Intuitive Middle Ear Testing

Our touchTymp impedance line is designed to turn tympanometry into an intuitive, efficient and enjoyable procedure. Improve your daily workflow with our full 10.4" touchscreen and a user-friendly interface that allows an easy change of parameters.





*„I can operate all features
within 3 clicks.*

touchTymp really is amazingly intuitive“

Dr. Michel Bloch, Cannes, France



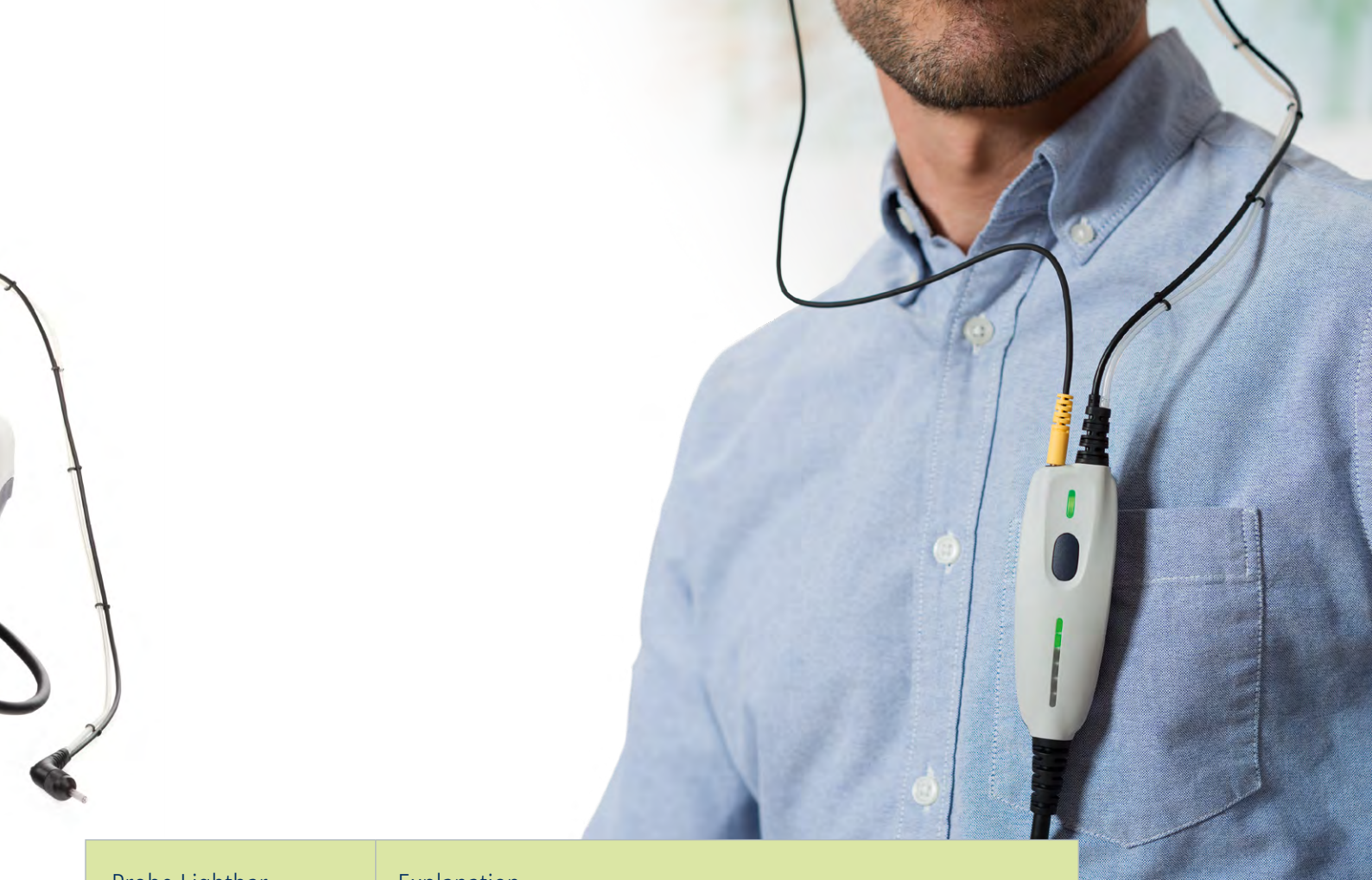

Probe Lightbar	Explanation	
	Tympanometry & Acoustic Reflex	Shows result: No Response
	Tympanometry & Acoustic Reflex	Shows result: Pass
	Acoustic Reflex	Stimulus is being given (additionally the last result is shown)
	Tympanometry	Lights up (rolling up) dependent on the values (normative box)

Focus on Your Patient

Experience full control of the probe while concentrating on your patient: touchTymp's easy-to-handle probes feature unique light bars to provide a real-time progression of the test. To exactly match your needs, the different touchTymp versions come with distinctive probes. touchTymp MI 24 and touchTymp MI 26 offer an ergonomic pen probe for screening tests. touchTymp MI 34 and touchTymp MI 36 feature a lightweight shoulder box probe for both screening and diagnostic purposes. It can easily be clipped to your patient's top for controlled handling.




With Audiometry for a Small Footprint

Our touchTymp MI 26 and touchTymp MI 36 both combine intuitive tympanometry and audiometry in one exceptional device. Their small footprint makes them perfect for limited desk space, and the ability to switch easily and smoothly between tests makes your workflows much more efficient and comfortable.



Choose between the touchTymp MI 26 for middle ear and audiometric screenings or the touchTymp MI 36 for added diagnostic capabilities. In addition to air conduction audiometry, the touchTymp MI 36 also offers bone conduction, which is available as an optional function for the touchTymp MI 26.

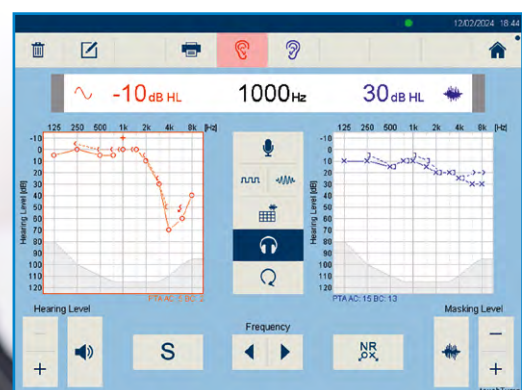

Tympanometry for Children has Never been so Exciting

We make immittance testing easy and fast by focusing the child's attention on an exciting animated car race for the duration of the tympanometry or reflex screening procedure. The colorful animation engages the child, reducing the likelihood of breaking the seal between probe tip and ear and having to retest. This saves you time and effort, while also making the screening procedure more enjoyable for your little patient!



The car race starts as soon as the test begins. Once the measurement is completed, the car reaches the finish line and the patient has won the race.



Comprehensive Test Protocols

Each touchTymp provides comprehensive standard protocols for immediate operation of screening and diagnostic tests:

Tests	touchTymp			
	MI 24	MI 26	MI 34	MI 36
Tympanometry				
226 Hz	✓	✓	✓	✓
678 Hz, 800 Hz			✓	✓
1000 Hz	optional upgrade	optional upgrade	optional upgrade	optional upgrade
RaceCar Animation 	✓	✓	✓	✓
Acoustic Reflexes				
Fixed (Screening)	✓	✓	✓	✓
Automatic (Threshold)	✓	✓	✓	✓
Ipsilateral	✓	✓	✓	✓
Contralateral	optional upgrade	optional upgrade	✓	✓
Pure tone	✓	✓	✓	✓
Noise			✓	✓
RaceCar Animation 	✓	✓	✓	✓
Reflex Decay				
Ipsilateral / Contralateral			✓	✓
ETF				
ETF Intact / Perforated			✓	✓
Audiometry				
Pure tone + warble tone		✓		✓
Air conduction		✓		✓
Bone conduction		optional upgrade		✓
Masking		✓		✓



Sessions PC Software

The touchTympanometer is fully supported by Sessions, so that test results can be transferred and managed on your PC. For our touchTympanometer, Sessions automatically synchronizes measurements ensuring a seamless transfer of your test data. With Sessions you can view, archive, manage the measurement data and print reports. The clean icon-based user interface enables an efficient workflow and high-speed data management. You can link up Sessions to your database with full compatibility or run it standalone. Alternatively, you can integrate it easily to your chosen EMR/ EHR solution via interfaces. This gives you the ability to create detailed reports that can be easily filed or printed. You can also create a "paperless" office by saving test results as PDFs for electronic filing or email.



Middle Ear Testing for All Ages

You and your patients are diverse – and so is our touchTympanometer. touchTympanometer is designed for the needs of all ages. Perform screening and diagnostic tests from infants patients to older adults.



Printing Made Easy

Save time by using the touchTympanometer's built-in printer: Just place the probe into its holder to instantly and automatically print test results and reports.




touchTymp – Components, Accessories & Software



touchTymp MI 24,
touchTymp MI 26
with Pen Probe




touchTymp MI 34,
touchTymp MI 36
with Shoulder Box
Probe

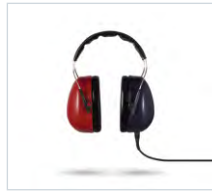

AC headphones
DD65 v2



BC headphone B71

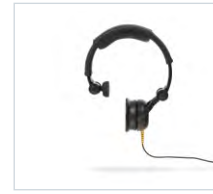

Contra headphone
DD45C




Eartip box



AC headphones
DD45



Talk back micro-
phone

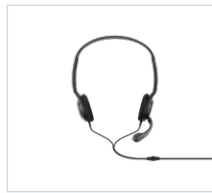

Monitor phone &
microphone headset




Sessions
PC software

Standard configuration vary according to touchTymp version choosen.

Disposables


Sanibel™ Supply is the exclusive supplier of MAICO touchTymp disposables. Use only Sanibel™ disposables to achive optimal test results.

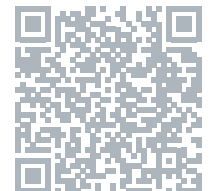

Watch our
touchTymp
videos



MAICO Diagnostics GmbH
Sickingenstr. 70 -71 · 10553 Berlin · Germany
Tel.: +49 30 / 70 71 46-50
sales@maico.biz · www.maico.biz

 Find your local distributor.



07/2025