



DR. BRILL + DR. STEINMANN
INSTITUT FÜR HYGIENE UND MIKROBIOLOGIE



05.06.2015

Test report C15L0096aPo

Evaluation of the effectiveness of Chemides Pulver

Test virus: poliovirus type 1 strain LSc-2ab

Method: EN 14476:2013/FprA1 2015

quantitative suspension test for the evaluation
of virucidal activity of chemical disinfectants and
antiseptics used in human medicine

Sponsor:

Chemi-Pharm AS
Pollu 132
EST – Tallinn 10917

Norderoog 2, D - 28259 Bremen
Tel.: +49 421-27819102, Fax: +49 421-2760283
info@brillhygiene.com, <http://www.brillhygiene.com>

TEST REPORT

1. Identification of test laboratory

Dr. Brill + Partner GmbH Institute for Hygiene and Microbiology, Norderoog 2, DE - 28259 Bremen

2. Identification of sample

Manufacturer	Chemi-Pharm AS
Name of product	Chemides Pulver
Product diluent recommended by the manufacturer	-
Batch number	30180215
Application	disinfectant for instruments and surfaces
Production date	18.02.2015
Expiry date	2017/02
Active compound (s) (kg)	peracetic acid in reaction mixture
Appearance, odour	white powder product specific
pH-values (in WSH)	1.75 %: 8.21 (20 °C)
Storage conditions	room temperature in the dark (area with restricted access)
Date of arrival in the laboratory	23.02.2015

3. Materials

3.1 Culture medium and reagents

- Dulbecco's Modified Eagle's Medium (DMEM, Biozym Scientific GmbH, catalogue no. 880021)
- fetal calf serum (Biochrom AG, article no. S 0115)
- 1.4 % formaldehyde solution (dilution of Roti®-Histofix 4 %, Carl Roth GmbH)
- Aqua bidest. (Fresenius Kabi Deutschland, article no. P2N 1636071)
- PBS (Invitrogen, article no. 18912-014)
- BSA (Sigma-Aldrich-Chemie GmbH, article no. CA-2153)
- sheep erythrocytes (Fiebig-Nährstofftechnik).

* Test procedure accredited according to DIN EN ISO/IEC 17025. Test report issued by Dr. Brill + Partner GmbH, Norderoog 2, DE – 28259 Bremen, Germany, Telephone +49. 421. 27819102, Telefax +49. 421. 2760283, www.brillhygiene.com. No copying or transmission, in whole or in part, of this test report without the explicit prior written permission. The test results exclusively apply to the tested samples. Information on measurement uncertainty on request. © Dr. Brill + Partner GmbH 2015

3.2 Virus and cells

The poliovirus type 1 strain LSc-2ab (Chiron-Behring) was obtained from PD Dr. Olaf Thraenhart, Eurovir, D-14943 Luckenwalde.

BGM cells (buffalo green monkey = permanent monkey kidney cell line; supplied by Prof. Dr. Lindl, Institut für angewandte Zellkultur, D-81669 München, Germany) were cultivated in a 175 cm² flask with Dulbecco's Modified Eagle's Medium (DMEM) and 10 % fetal calf serum (FCS).

The cells (passage 38) were inspected regularly for morphological alterations and for contamination by mycoplasmas. No morphological alterations of cells and no contamination by mycoplasmas could be detected.

3.3 Apparatus, glassware and small items of equipment

- CO₂ incubator, Nunc GmbH & Co. KG, model QWJ 350
- Agitator (Vortex Genie Mixer, type G 560E)
- pH measurement 315i (WTW, article no. 2A10-100)
- Centrifuge (Sigma-Aldrich-Chemie GmbH, type 113)
- Microscope (Olympus, type CK 30)
- Centrifuge 5804 R (Eppendorf AG)
- Water bath (JULABO, Julabo U 3)
- Adjustable and fixed-volume pipettes (Eppendorf AG)
- Polystyrol 96-well microtitre plate (Nunc GmbH & Co. KG, Wiesbaden)
- Cell culture flask (Nunc GmbH & Co. KG, Wiesbaden)
- Sealed test tubes (Sarstedt AG & Co., Nümbrecht).

* Test procedure accredited according to DIN EN ISO/IEC 17025. Test report issued by Dr. Brill + Partner GmbH, Norderoog 2, DE – 28259 Bremen, Germany, Telephone +49. 421. 27819102, Telefax +49. 421. 2760283, www.brillhygiene.com. No copying or transmission, in whole or in part, of this test report without the explicit prior written permission. The test results exclusively apply to the tested samples. Information on measurement uncertainty on request. © Dr. Brill + Partner GmbH 2015

4. Experimental conditions

Test temperature	20 °C ± 1.0 °C
Concentration of test product	1.75 %, 1.0 %, 0.5 % and 0.1 % (demonstration of non-active range) solutions
Appearance of product dilutions	no precipitation
Contact times	10, 15, 30 and 60 minutes
Interfering substance	3.0 g/l bovine serum albumin + 3.0 g/l erythrocytes (dirty conditions EN 14476:2013)
Procedure to stop action of disinfectant	immediate dilution
Diluent	water of standardised hardness (WSH)
Stability of product in the mix with virus and interfering substance (1.75 % solution)	no flocculation, no precipitation
Virus strain	poliovirus type 1 strain LSc-2ab (Chiron-Behring)
Date of testing	23.02.2015 – 05.06.2015
End of testing	05.06.2015

5. Methods

5.1 Preparation of test virus suspension

For preparation of test virus suspension according to EN 5.4.1 *BGM cells* were infected with a multiplicity of infection of 0.1 at 37 °C. After cells showed a cytopathic effect, they were subjected to a threefold freeze/thaw procedure followed by a low speed centrifugation in order to sediment cell debris. After aliquotation of the supernatant, test virus suspension was stored at -80 °C.

5.2 Preparation of disinfectant (dilutions)

A working solution of the test product was prepared by dissolving 17.5 g of the powder (taking into account a factor of 1.25 due to addition of 1 part test virus and 1 part interfering substance to 8 parts of disinfectant) in pre-warmed WSH. Further dilutions (1.0 %, 0.5 % and 0.1 %) were prepared out of this working solution immediately before the inactivation tests.

* Test procedure accredited according to DIN EN ISO/IEC 17025. Test report issued by Dr. Brill + Partner GmbH, Norderoog 2, DE – 28259 Bremen, Germany, Telephone +49. 421. 27819102, Telefax +49. 421. 2760283, www.brillhygiene.com. No copying or transmission, in whole or in part, of this test report without the explicit prior written permission. The test results exclusively apply to the tested samples. Information on measurement uncertainty on request. © Dr. Brill + Partner GmbH 2015

5.3 Infectivity assay

Infectivity was determined as endpoint titration according to EN 5.5 transferring 0.1 ml of each dilution into eight wells of a microtitre plate to 0.1 ml of freshly trypsinised *BGM cells* ($10\text{--}15 \times 10^3$ cells per well), beginning with the highest dilution. Microtitre plates were incubated at 37 °C in a 5 % CO₂-atmosphere. The cytopathic effect was read by using an inverted microscope after seven days. Calculation of the infective dose TCID₅₀/ml was calculated with the method of Spearman (2) and Kärber (3) with the following formula:

$$-\log_{10}\text{TCID}_{50} = X_0 - 0.5 + \sum r/n$$

meaning

X_0 = $\log_{10}$ of the lowest dilution with 100 % positive reaction

r = number of pos. determinations of lowest dilution step with 100 % positive and all higher positive dilution steps

n = number of determinations for each dilution step.

5.4 Calculation and verification of virucidal activity

The virucidal activity of the test disinfectant was evaluated by calculating the decrease in titre in comparison with the control titration without disinfectant. The difference is given as reduction factor (RF).

According to the EN 14476:2013, a disinfectant or a disinfectant solution at a particular concentration is having virus-inactivating efficacy if the titre is reduced at least by four $\log_{10}$ steps within the recommended exposure period. This corresponds to an inactivation of ≥ 99.99 %.

5.5 Inactivation assay

Determination of virucidal activity has been carried out in accordance to EN 5.5. The test product was examined as 1.75 %, 1.0 %, 0.5 % and 0.1 % (demonstration of non-active range) solutions in WSH at 20 °C according to EN 14476:2013. 10, 15, 30 and 60 minutes were chosen as contact times.

Immediately at the end of a chosen contact time, activity of the disinfectant was stopped by dilution to 10^{-8} .

Titration of the virus control were performed after the longest exposure time (EN 5.5.7).

* Test procedure accredited according to DIN EN ISO/IEC 17025. Test report issued by Dr. Brill + Partner GmbH, Norderoog 2, DE – 28259 Bremen, Germany, Telephone +49. 421. 27819102, Telefax +49. 421. 2760283, www.brillhygiene.com. No copying or transmission, in whole or in part, of this test report without the explicit prior written permission. The test results exclusively apply to the tested samples. Information on measurement uncertainty on request. © Dr. Brill + Partner GmbH 2015

Furthermore, a cell control (only addition of medium) was incorporated.

Inactivation tests were carried out in sealed test tubes in a water bath at $20\text{ °C} \pm 1.0\text{ °C}$. Aliquots were retained after appropriate exposure times and residual infectivity was determined.

5.6 Determination of cytotoxicity

Determination of cytotoxicity was performed according to EN 5.5.4.1.

5.7 Cell sensitivity to virus

For the control of cell sensitivity to virus two parts by volume hard water were mixed with eight parts by volume of the lowest apparently non-cytotoxic dilution of the product. This mixture or PBS as control was added to a volume of double concentrated cell suspension. After 1 h at 37 °C the cells were centrifuged and re-suspended in cell culture medium (EN 5.5.4.2b).

Finally, a comparative titration of the test virus suspension was performed on the pre-treated (disinfectant) and non-pre-treated (PBS) cells as described above.

5.8 Control of efficacy for suppression of disinfectant's activity

Furthermore, a control of efficiency for suppression of disinfectant's activity was included (EN 5.5.5).

5.9 Reference virus inactivation test

As reference for test validation a 0.7 % formaldehyde solution according to EN 5.5.6 was included. 5, 15, 30 and 60 minutes were chosen as contact times. In addition, cytotoxicity of formaldehyde test solution was determined following EN 5.5.6.2 with dilutions up to 10^{-5} .

* Test procedure accredited according to DIN EN ISO/IEC 17025. Test report issued by Dr. Brill + Partner GmbH, Norderoog 2, DE – 28259 Bremen, Germany, Telephone +49. 421. 27819102, Telefax +49. 421. 2760283, www.brillhygiene.com. No copying or transmission, in whole or in part, of this test report without the explicit prior written permission. The test results exclusively apply to the tested samples. Information on measurement uncertainty on request. © Dr. Brill + Partner GmbH 2015

6. Verification of the methodology

The following criteria as mentioned in EN 5.7 were fulfilled:

- a) The titre of the test virus suspension allowed the determination of a $\geq 4 \log_{10}$ reduction (maximal virus reduction $\geq 5.63 \pm 0.31$).
- b) The difference of the logarithmic titre of the virus control minus the logarithmic titre of the test virus in the reference inactivation test (see 6.6.7) was 1.63 ± 0.62 (between 0.5 - 2.5) after 30 min and 2.75 ± 0.50 (between 2.0 - 4.5) after 60 min for poliovirus type 1.
- c) The test product (1.75 % and 1.0 % solutions) showed cytotoxicity in the 1:10 dilutions thus allowing the detection of a $4 \log_{10}$ reduction of virus titre.
- d) The comparative titration on pre-treated (disinfectant) and non-pre-treated (PBS) *BGM cells* showed no significant difference ($< 1 \log_{10}$; EN 5.7) of virus titre: 7.63 ± 0.25 (PBS) versus 8.00 ± 0.38 (1:100 dilutions of disinfectant, 1.75 %) $\log_{10}\text{TCID}_{50}/\text{ml}$.
- e) The control of efficacy for suppression of disinfectant's activity (1.75 %) showed no decrease ($< 0.5 \log_{10}$; EN 5.5.5.1) in virus titre (7.63 ± 0.25 versus $7.75 \pm 0.35 \log_{10}\text{TCID}_{50}/\text{ml}$).
- f) One concentration demonstrated a $4 \log_{10}$ reduction and (at least) one concentration demonstrated a $\log_{10}$ reduction of less than 4.

Since all criteria according EN 5.7 were fulfilled, examination with poliovirus type 1 according to EN 14476:2013 is valid.

* Test procedure accredited according to DIN EN ISO/IEC 17025. Test report issued by Dr. Brill + Partner GmbH, Norderoog 2, DE – 28259 Bremen, Germany, Telephone +49. 421. 27819102, Telefax +49. 421. 2760283, www.brillhygiene.com. No copying or transmission, in whole or in part, of this test report without the explicit prior written permission. The test results exclusively apply to the tested samples. Information on measurement uncertainty on request. © Dr. Brill + Partner GmbH 2015



7. Results

Results of examination are shown in tables 1 to 9. Tables 1 to 8 demonstrate the raw data, whereas table 9 (a+b) gives a summary of results.

The test product as 1.75 % solution was able to inactivate poliovirus type 1 after 15 minutes in this quantitative suspension test (Tables 1 and 2). The reduction factors were $\geq 5.00 \pm 0.58$ and $\geq 4.50 \pm 0.57$ at this time point. The mean value was $\geq 4.75 \pm 0.41$. This corresponded to an inactivation of ≥ 99.99 %.

The test product as 1.0 % solution was able to inactivate poliovirus type 1 after 30 minutes in this quantitative suspension test (Table 3). The reduction factor was $\geq 4.25 \pm 0.63$.

Tested as 0.5 % solution, the test product was not sufficient active within 60 minutes of exposure time (Table 4).

In addition, tested as 0.1 % solution, the test product was not active within 30 minutes of exposure time (Table 5).

8. Conclusion

The instrument disinfectant Chemides Pulver tested as 1.75 % solution demonstrated effectiveness against poliovirus type 1 after an exposure time of 15 minutes under dirty conditions. The 1.0 % solution was active within 30 minutes. Therefore, the instrument disinfectant Chemides Pulver can be declared as active against poliovirus type 1 as follows:

1.75 %	15 minutes
1.0 %	30 minutes

Bremen, 05.06.2015

- Dr. Jochen Steinmann -
Scientific Director



* Test procedure accredited according to DIN EN ISO/IEC 17025. Test report issued by Dr. Brill + Partner GmbH, Norderoog 2, DE - 28259 Bremen, Germany, Telephone +49. 421. 27819102, Telefax +49. 421. 2760283, www.brillhygiene.com. No copying or transmission, in whole or in part, of this test report without the explicit prior written permission. The test results exclusively apply to the tested samples. Information on measurement uncertainty on request. © Dr. Brill + Partner GmbH 2015



9. Quality control

The Quality Assurance of the results was maintained by performing the determination of the virus-inactivating properties of the disinfectant in accordance with Good Laboratory Practice regulations:

- 1) Chemicals Act of Germany, Appendix 1, dating of 01.08 1994 (BGBl. I, 1994, page 1703). Appendix revised at 14. 05. 1997 (BGBl. I, 1997, page 1060).
- 2) OECD Principles of Good Laboratory Practice (revised 1997); OECD Environmental Health and Safety Publications; Series on Principles of Good Laboratory Practice and Compliance Monitoring – Number 1. Environment Directorate, Organization for Economic Co-operation and Development, Paris 1998.

The plausibility of the results was additionally confirmed by controls incorporated in the inactivation assays.

10. Records to be maintained

All testing data, protocol, protocol modifications, the final report, and correspondence between Dr. Brill + Partner GmbH and the sponsor will be stored in the archives at Dr. Brill + Partner GmbH.

The use of the Dr. Brill + Partner GmbH name, logo or any other representation of Dr. Brill + Partner GmbH, other than distribution of this report in it's entirely, without the written approval of Dr. Brill + Partner GmbH is prohibited. In addition, Dr. Brill + Partner GmbH may not be referred to in any form of promotional materials, press releases, advertising or similar materials (whether by print, broadcast, communication or electronic means) without the express permission of Dr. Brill + Partner GmbH.

The test results in this test report relate only to the items examined.

* Test procedure accredited according to DIN EN ISO/IEC 17025. Test report issued by Dr. Brill + Partner GmbH, Norderoog 2, DE – 28259 Bremen, Germany, Telephone +49. 421. 27819102, Telefax +49. 421. 2760283, www.brillhygiene.com. No copying or transmission, in whole or in part, of this test report without the explicit prior written permission. The test results exclusively apply to the tested samples. Information on measurement uncertainty on request. © Dr. Brill + Partner GmbH 2015



11. Literature

1. EN 14476:2013/FprA1 2015: Chemical disinfectants and antiseptics – Quantitative suspension test for the evaluation of virucidal activity of chemicals disinfectants and antiseptics in human medicine test - Test method and requirements (phase 2, step 1)
2. Spearman, C.: The method of 'right or wrong cases' (constant stimuli) without Gauss's formulae.
Brit J Psychol; 2 1908, 227-242
3. Kärber, G.: Beitrag zur kollektiven Behandlung pharmakologischer Reihenversuche. [A contribution to the collective treatment of a pharmacological experimental series].
Arch Exp Path Pharmac; 162, 1931, 480-487

* Test procedure accredited according to DIN EN ISO/IEC 17025. Test report issued by Dr. Brill + Partner GmbH, Norderoog 2, DE – 28259 Bremen, Germany, Telephone +49. 421. 27819102, Telefax +49. 421. 2760283, www.brillhygiene.com. No copying or transmission, in whole or in part, of this test report without the explicit prior written permission. The test results exclusively apply to the tested samples. Information on measurement uncertainty on request. © Dr. Brill + Partner GmbH 2015



Appendix:

Legend to the Tables

Table 1:	Raw data for Chemides Pulver (1.75 %) tested against poliovirus type 1 (1 st assay)
Table 2:	Raw data for Chemides Pulver (1.75 %) tested against poliovirus type 1 (2 nd assay)
Table 3:	Raw data for Chemides Pulver (1.0 %) tested against poliovirus type 1
Table 4:	Raw data for Chemides Pulver (0.5%) tested against poliovirus type 1
Table 5:	Raw data for Chemides Pulver (0.1 %) tested against poliovirus type 1
Table 6:	Raw data for formaldehyde solution (0.7 %) tested against poliovirus type 1
Table 7:	Raw data for control of efficacy for suppression of disinfectant's activity (1.75 %)
Table 8:	Raw data (poliovirus type 1) for cell sensitivity (1.75 %)
Table 9 (a+b):	Summary of results with Chemides Pulver and poliovirus type 1

Legend to the Figures

Figure 1:	Virus-inactivating properties of Chemides Pulver (1.75 %)
Figure 2:	Virus-inactivating properties of formaldehyde (0.7 %)

* Test procedure accredited according to DIN EN ISO/IEC 17025. Test report issued by Dr. Brill + Partner GmbH, Norderoog 2, DE – 28259 Bremen, Germany, Telephone +49. 421. 27819102, Telefax +49. 421. 2760283, www.brillhygiene.com. No copying or transmission, in whole or in part, of this test report without the explicit prior written permission. The test results exclusively apply to the tested samples. Information on measurement uncertainty on request. © Dr. Brill + Partner GmbH 2015



Table 1: Raw data for Chemides Pulver (1.75 %) tested against poliovirus type 1 at 20 °C (quantal test; 8 wells) (3843) (1st assay)

Product	Concentration	Interfering substance	Contact time (min)	Dilutions (log ₁₀)								
				1	2	3	4	5	6	7	8	9
test product	1.75%	dirty conditions	10	n.d.	n.d.	n.d.	n.d.	n.d.	n.d.	n.d.	n.d.	n.d.
			15	tttt tttt	4004 0444	0000 0000	0000 0000	0000 0000	0000 0000	0000 0000	n.d.	n.d.
			30	tttt tttt	0000 0000	0000 0000	0000 0000	0000 0000	0000 0000	0000 0000	n.d.	n.d.
			60	n.d.	n.d.	n.d.	n.d.	n.d.	n.d.	n.d.	n.d.	n.d.
test product cytotoxicity	1.75%	dirty conditions	n.a.	tttt tttt	0000 0000	0000 0000	0000 0000	0000 0000	n.d.	n.d.	n.d.	n.d.
virus control	n.a.	dirty conditions	0	n.d.	n.d.	n.d.	n.d.	n.d.	n.d.	n.d.	n.d.	n.d.
			60	4444 4444	4444 4444	4444 4444	4444 4444	4444 4444	4444 4444	0404 4040	4000 0000	0000 0000

n.a. = not applicable
n.d. = not done

0 = no virus present; t = cytotoxic
1 to 4 = virus present (degree of CPE in 8 cell culture units) (wells of microtitre plates)



Table 2: Raw data for Chemides Pulver (1.75 %) tested against poliovirus type 1 at 20 °C (quantal test; 8 wells) (3925) (2nd assay)

Product	Concentration	Interfering substance	Contact time (min)	Dilutions (log ₁₀)								
				1	2	3	4	5	6	7	8	9
test product	1.75%	dirty conditions	10	tttt	4444	4444	4444	0000	0000	0000	n.d.	n.d.
				tttt	4444	4444	4444	0040	0000	0000	n.d.	n.d.
			15	tttt	4004	0400	0000	0000	0000	0000	n.d.	n.d.
				tttt	4440	0000	0000	0000	0000	0000	n.d.	n.d.
test product cytotoxicity	1.75%	dirty conditions	30	n.d.	n.d.	n.d.	n.d.	n.d.	n.d.	n.d.	n.d.	
				n.d.	n.d.	n.d.	n.d.	n.d.	n.d.	n.d.	n.d.	
			60	n.d.	n.d.	n.d.	n.d.	n.d.	n.d.	n.d.	n.d.	n.d.
				n.d.	n.d.	n.d.	n.d.	n.d.	n.d.	n.d.	n.d.	n.d.
virus control	n.a.	dirty conditions	n.a.	tttt	0000	0000	0000	0000	n.d.	n.d.	n.d.	
				tttt	0000	0000	0000	0000	0000	0000	0000	0000
			0	4444	4444	4444	4444	4444	4444	0003	0000	0000
				4444	4444	4444	4444	4444	4444	3000	0030	0000
			60	4444	4444	4444	4444	4444	4444	0003	0000	0000
				4444	4444	4444	4444	4444	4444	0000	0030	0000

n.a. = not applicable
n.d. = not done

0 = no virus present; t = cytotoxic

1 to 4 = virus present (degree of CPE in 8 cell culture units) (wells of microtitre plates)



Table 3: Raw data for Chemides Pulver (1.0 %) tested against poliovirus type 1 at 20 °C (quantal test; 8 wells) (3925)

Product	Concentration	Interfering substance	Contact time (min)	Dilutions (log ₁₀)									
				1	2	3	4	5	6	7	8	9	
test product	1.0%	dirty conditions	10	n.d.	n.d.	n.d.	n.d.	n.d.	n.d.	n.d.	n.d.	n.d.	
			15	tttt	4444	4444	4444	4344	4400	0000	0000	n.d.	
			30	tttt	4444	4444	0000	0000	0000	0000	0000	n.d.	
			60	tttt	0440	0000	0000	0000	0000	0000	0000	n.d.	
test product cytotoxicity	1.0%	dirty conditions	n.a.	tttt	0000	0000	0000	0000	0000	n.d.	n.d.	n.d.	
virus control	n.a.	dirty conditions	0	4444	4444	4444	4444	4444	4444	4444	0003	0000	0000
			60	4444	4444	4444	4444	4444	4444	4444	0003	0000	0030

n.a. = not applicable
n.d. = not done

0 = no virus present; t = cytotoxic
1 to 4 = virus present (degree of CPE in 8 cell culture units) (wells of microtitre plates)



Table 4: Raw data for Chemides Pulver (0.5 %) tested against poliovirus type 1 at 20 °C (quantal test; 8 wells) (3925)

Product	Concentration	Interfering substance	Contact time (min)	Dilutions (log ₁₀)									
				1	2	3	4	5	6	7	8	9	
test product	0.5%	dirty conditions	10	n.d.	n.d.	n.d.	n.d.	n.d.	n.d.	n.d.	n.d.	n.d.	n.d.
			15	n.d.	n.d.	n.d.	n.d.	n.d.	n.d.	n.d.	n.d.	n.d.	n.d.
			30	4444 4444	4444 4444	4444 4444	4444 4444	4444 4444	0003 3404	0000 0400	0000 0000	n.d.	
			60	4444 4444	4444 4444	4444 4444	4430 4444	0000 0000	0000 0000	0000 0000	0000 0000	n.d.	
test product cytotoxicity	0.5%	dirty conditions	n.a.	0000 0000	0000 0000	0000 0000	0000 0000	n.d.	n.d.	n.d.	n.d.		
virus control	n.a.	dirty conditions	60	4444 4444	4444 4444	4444 4444	4444 4444	4444 4444	4444 4444	0003 0000	0000 0030	0000 0000	

n.a. = not applicable
n.d. = not done

0 = no virus present; t = cytotoxic

1 to 4 = virus present (degree of CPE in 8 cell culture units) (wells of microtitre plates)



Table 5: Raw data for Chemides Pulver (0.1 %) tested against poliovirus type 1 at 20 °C (quantal test; 8 wells) (3925)

Product	Concentration	Interfering substance	Contact time (min)	Dilutions (log ₁₀)								
				1	2	3	4	5	6	7	8	9
test product	0.1%	dirty conditions	10	n.d.	n.d.	n.d.	n.d.	n.d.	n.d.	n.d.	n.d.	n.d.
			15	n.d.	n.d.	n.d.	n.d.	n.d.	n.d.	n.d.	n.d.	n.d.
			30	4444 4444	4444 4444	4444 4444	4344 4244	4444 4444	4444 4444	0400 0000	0000 0000	n.d. n.d.
			60	n.d.	n.d.	n.d.	n.d.	n.d.	n.d.	n.d.	n.d.	n.d.
test product cytotoxicity	0.1%	dirty conditions	n.a.	0000 0000	0000 0000	0000 0000	0000 0000	0000 0000	n.d.	n.d.	n.d.	
virus control	n.a.	dirty conditions	60	4444 4444	4444 4444	4444 4444	4444 4444	4444 4444	4444 4444	0003 0000	0000 0030	0000 0000

n.a. = not applicable
n.d. = not done

0 = no virus present; t = cytotoxic

1 to 4 = virus present (degree of CPE in 8 cell culture units) (wells of microtitre plates)



Table 6: Raw data for formaldehyde solution (0.7 %) tested against poliovirus type 1 at 20 °C (quantal test; 8 wells) (3925)

Product	Concentration	Interfering substance	Contact time (min)	Dilutions (log ₁₀)								
				1	2	3	4	5	6	7	8	9
formaldehyde	0.7% (m/V)	PBS	5	tttt	tttt	4444	4444	4444	4404	0000	0000	n.d.
			15	tttt	tttt	4444	4444	4444	4400	0000	0000	n.d.
			30	tttt	tttt	4444	4444	4440	0000	0000	0000	n.d.
			60	tttt	tttt	4444	4304	0000	0000	0000	0000	n.d.
formaldehyde cytotoxicity	0.7% (m/V)	PBS	n.a.	tttt	tttt	0000	0000	0000	n.d.	n.d.	n.d.	n.d.
virus control	n.a.	PBS	0	n.d.	n.d.	n.d.	n.d.	n.d.	n.d.	n.d.	n.d.	n.d.
			60	4444	4444	4444	4444	4444	4444	0040	0000	0000

n.a. = not applicable
n.d. = not done

0 = no virus present; t = cytotoxic
1 to 4 = virus present (degree of CPE in 8 cell culture units) (wells of microtitre plates)



Deutsche
Akkreditierungsstelle
D-PL 13127-01-02



Zustimmung der Länder
für die Akkreditierung
bei Anträgen und
Beschwerden
ZLG-AP-306.10.31



Table 7: Raw data for control of efficacy for suppression of disinfectant's activity (1.75 %) (3925)

Product	Interfering substance	dilutions (log ₁₀)								
		1	2	3	4	5	6	7	8	9
test product	PBS	n.d.	n.d.	n.d.	n.d.	n.d.	n.d.	n.d.	n.d.	n.d.
test product	clean conditions	n.d.	n.d.	n.d.	n.d.	n.d.	n.d.	n.d.	n.d.	n.d.
test product	dirty conditions	tttt tttt	4444 4444	4444 4444	4444 4444	4444 4444	4444 4444	4000 0000	0000 0000	n.d.

n.a. = not applicable
n.d. = not done

0 = no virus present; t = cytotoxic

1 to 4 = virus present (degree of CPE in 8 cell culture units) (wells of microtitre plates)



Table 8: Raw data (poliovirus type 1) for cell sensitivity (1.75 %) (3925)

Product	Dilution	Dilutions (log ₁₀)								
		1	2	3	4	5	6	7	8	9
PBS	-	4444 4444	4444 4444	4444 4444	4444 4444	4444 4444	4444 4444	4000 0000	0000 0000	n.d.
test product	1:10	n.d.	n.d.	n.d.	n.d.	n.d.	n.d.	n.d.	n.d.	n.d.
test product	1:100	4444 4444	4444 4444	4444 4444	4444 4444	4444 4444	4444 4444	0044 4004	0000 0000	n.d.
test product	1:1,000	n.d.	n.d.	n.d.	n.d.	n.d.	n.d.	n.d.	n.d.	n.d.

n.a. = not applicable

0 = no virus present; t = cytotoxic

1 to 4 = virus present (degree of CPE in 8 cell culture units) (wells of microtitre plates)

n.d. = not done



Table 9a: Summary of results with Chemides Pulver and poliovirus type 1

Product	Con- centration	Interfering substance	Level of cytotoxicity	log ₁₀ TCID ₅₀ /ml aftermin					> 4 log ₁₀ reduction after ... min
				1	10	15	30	60	
test product	1.75%	dirty	2.50	n.d.	n.d.	≤3.13±0.37	≤2.50±0.00	n.d.	15 (RF ≥ 5.00±0.58)
test product	1.75%	dirty	2.50	n.d.	5.63±0.25	≤3.25±0.44	n.d.	n.d.	15 (RF ≥ 4.50±0.57)
test product	1.0%	dirty	2.50	n.d.	n.d.	6.75±0.33	3.50±0.52	≤2.50±0.00	30 (RF ≥ 4.25±0.63)
test product	0.5%	dirty	1.50	n.d.	n.d.	n.d.	7.13±0.45	5.38±0.25	> 60
test product	0.1%	dirty	1.50	n.d.	n.d.	n.d.	7.63±0.25	n.d.	> 30

n.a. = not applicable n.d. = not done




Table 9b: Summary of results with Chemides Pulver and poliovirus type 1

Product	Con- centration	Interfering substance	Level of cytotoxicity	log ₁₀ TCID ₅₀ /ml aftermin					> 4 log ₁₀ reduction after ... min
				0	5	15	30	60	
formaldehyde	0.7% (w/v)	PBS	3.50	n.d.	7.38±0.41	6.88±0.37	6.38±0.49	5.25±0.33	> 60
virus contr.	n.a.	PBS	n.a.	n.d.	n.d.	n.d.	n.d.	8.00±0.38	n.a.
virus contr. 1	n.a.	dirty	n.a.	n.d.	n.d.	n.d.	n.d.	8.13±0.45	n.a.
virus contr. 2	n.a.	dirty	n.a.	7.88±0.41	n.d.	n.d.	n.d.	7.75±0.35	n.a.
suppression control	1.75%	dirty	2.50	n.d.	n.d.	n.d.	7.63±0.25	n.d.	n.a.
sens.control PBS	n.a.	dirty	n.a.	n.d.	n.d.	n.d.	n.d.	7.63±0.25	n.a.
sens. control test product	1.75% → 1:100	dirty	n.a.	n.d.	n.d.	n.d.	n.d.	8.00±0.38	n.a.

n.a. = not applicable n.d. = not done sens. = sensitivity




Figure 1: Virus-inactivating properties of Chemides Pulver (1.75 %)

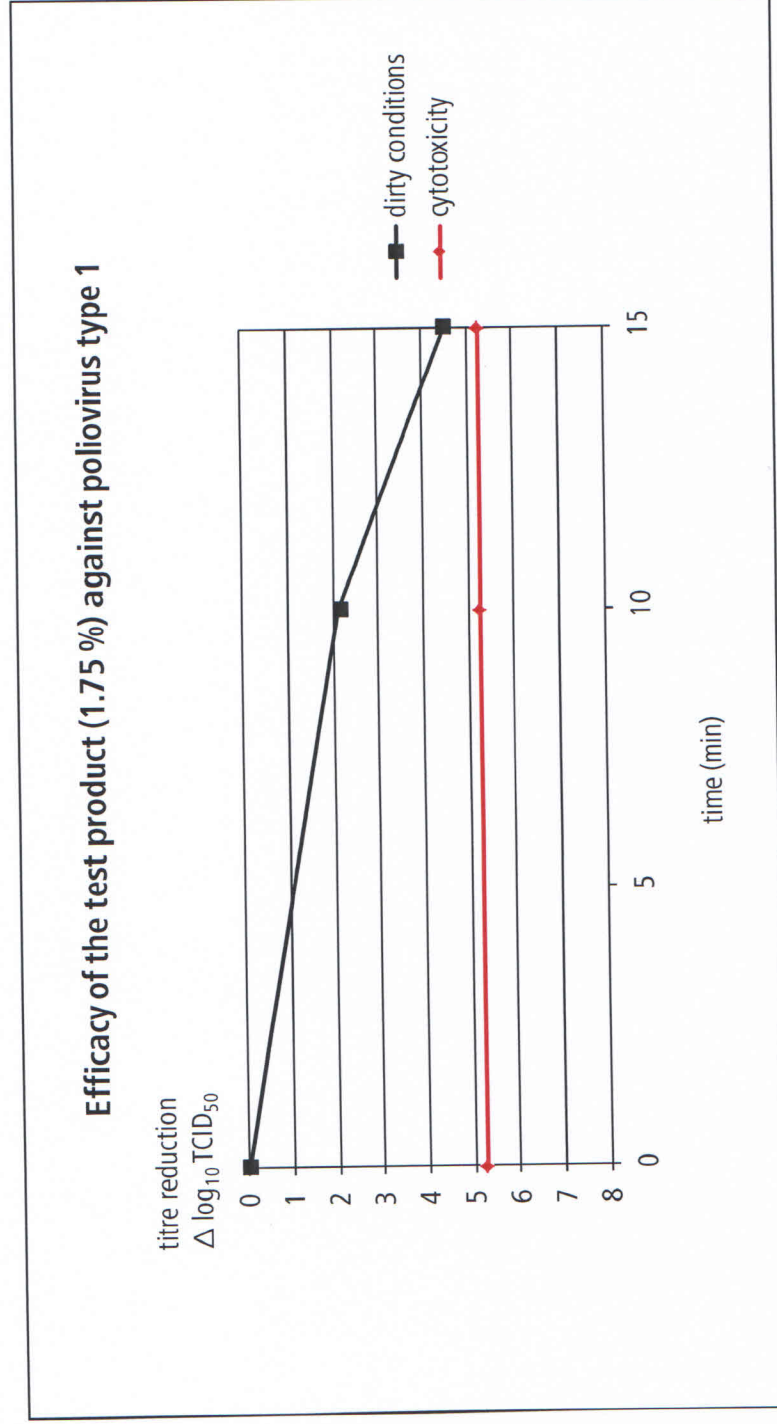

* Test procedure accredited according to DIN EN ISO/IEC 17025. Test report issued by Dr. Brill + Partner GmbH, Norderoog 2, DE – 28259 Bremen, Germany, Telephone +49. 421. 27819102, Telefax +49. 421. 2760283, www.brillhygiene.com. No copying or transmission, in whole or in part, of this test report without the explicit prior written permission. The test results exclusively apply to the tested samples. Information on measurement uncertainty on request. © Dr. Brill + Partner GmbH 2015



Figure 2: Virus-inactivating properties of formaldehyde (0.7 %) (titre reduction $\Delta \log_{10} \text{TCID}_{50}$)

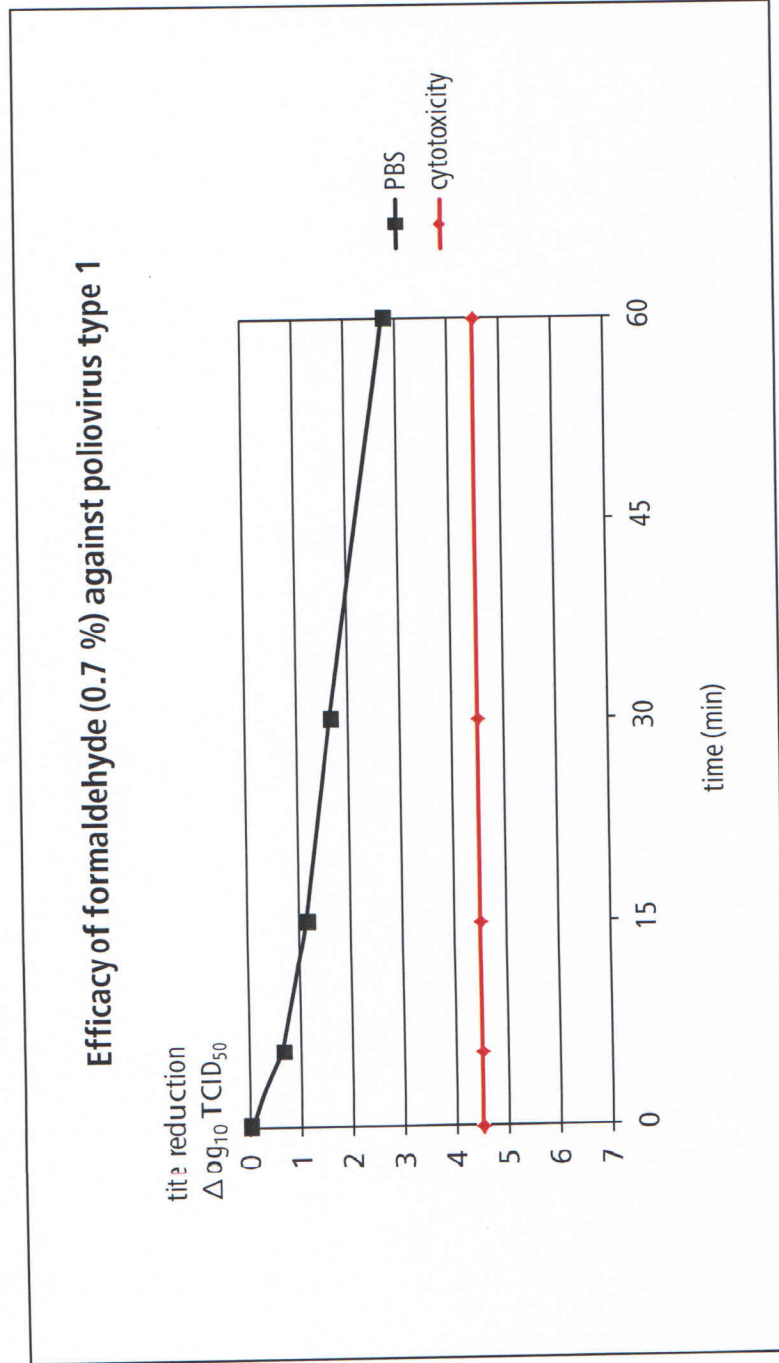

* Test procedure accredited according to DIN EN ISO/IEC 17025. Test report issued by Dr. Brill + Partner GmbH, Norderoog 2, DE – 28259 Bremen, Germany, Telephone +49. 421. 27819102, Telefax +49. 421. 2760283, www.brillhygiene.com. No copying or transmission, in whole or in part, of this test report without the explicit prior written permission. The test results exclusively apply to the tested samples. Information on measurement uncertainty on request. © Dr. Brill + Partner GmbH 2015

