



ASAP^{LP} ASAP[®]

ASPIRATION CATHETERS AND KITS

Your Options for
Effective and Quick
Clot Aspiration

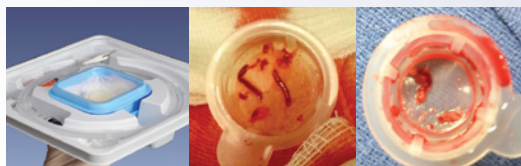

ASAPLP™ AND ASAP® ASPIRATION CATHETERS KEY FEATURES

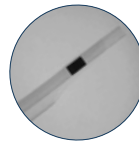

EFFECTIVE

- 100% stainless steel braided, kink resistant 145 cm low profile aspiration catheter
- Improved tip radiopacity due to radiopaque platinum iridium marker band
- Bright colored wire exchange lumen

QUICK ASPIRATION

- Fits into 6F guide catheter of ≥ 0.066 " with ease
- A 4F catheter designed for maximum aspiration
- Long 30 cm hydrophilic coating promoting deliverability
- Longer 20 cm rapid exchange lumen designed for better wire trackability

OPTIONS

- Packaged with stiffening stylet



EFFECTIVE

- Single extrusion, 100% stainless steel braided, kink resistant 140 cm aspiration catheter
- Tapered tip for smooth catheter transition

QUICK ASPIRATION

- Large aspiration lumen to facilitate quick aspiration
- 3 depth-positioning nonradiopaque markers at 90 cm, 100 cm, 110 cm
- All-in-one kit with all the necessary components to perform a catheter aspiration procedure

OPTIONS

- Standalone MiniStop® waste basin with lid and 70-micron filter



CONVENIENCE AT YOUR FINGERTIPS


Both catheters are available in the ASAP Kit with all the necessary components to perform a catheter aspiration procedure safely, quickly, and effectively. With our all-in-one tray, simply remove the lid to begin the procedure.

KIT COMPONENTS


- (1) Rapid exchange kink resistant ASAPLP or ASAP aspiration catheter
- (2) Merit VacLok® 30 mL aspiration syringes
- (1) 4 mL RXP® flushing syringe
- (2) Filter baskets, 70 micron
- (1) MicroStop™ closed waste basin, splatter proof, 60 mL fluid capacity
- (1) 21.5 cm tubing set with stopcock

A FAST, HANDY SOLUTION

The ASAP kits are designed for the quick removal of fresh, soft emboli and thrombi from vessels of the arterial system.


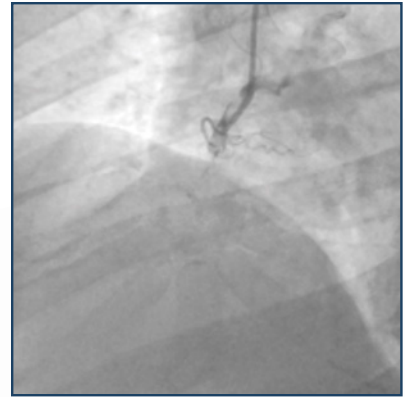


Pre Aspiration

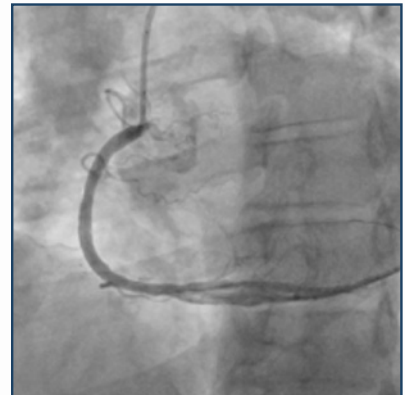

Post Aspiration

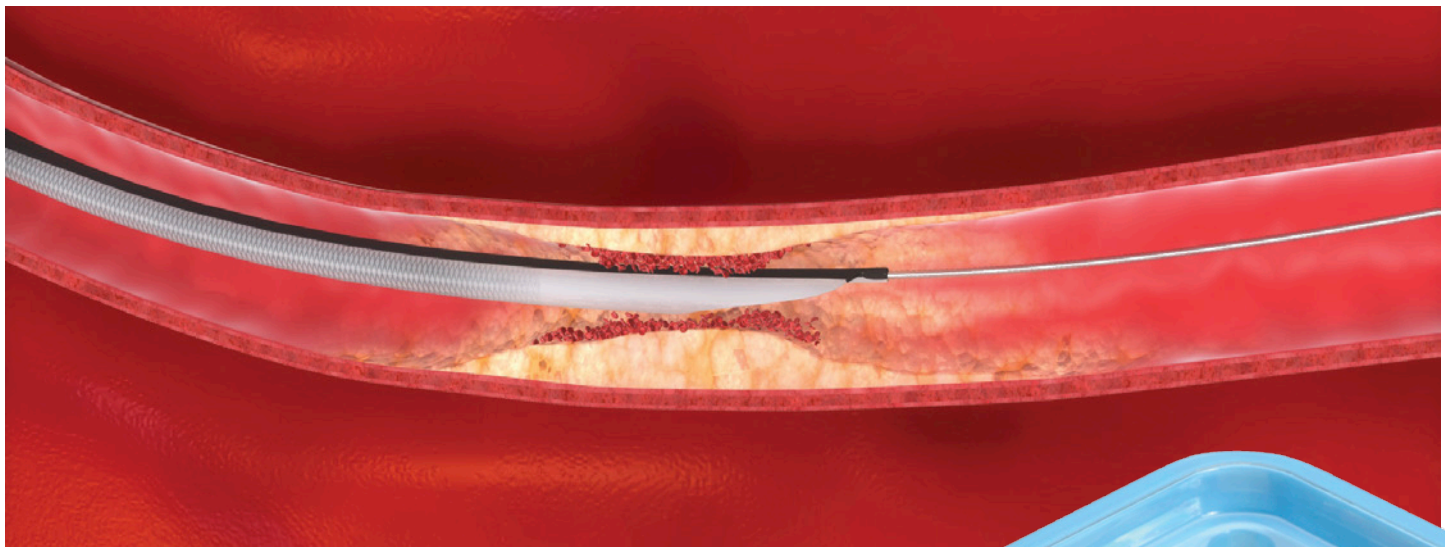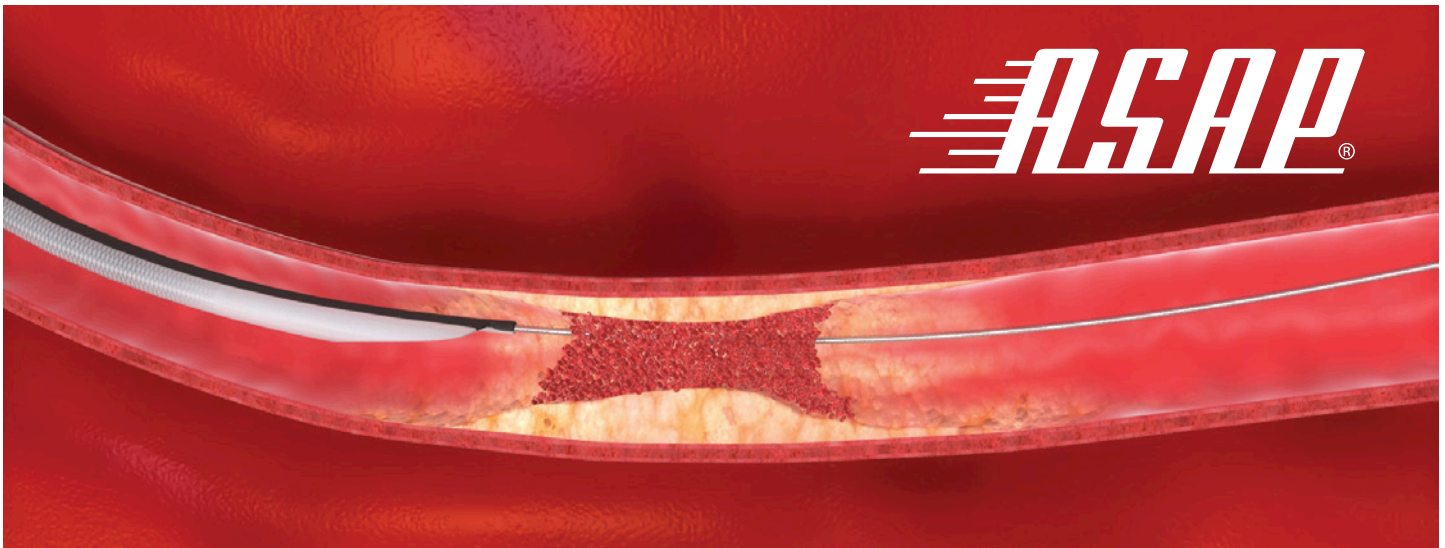

The ASAP logo is located in the top right corner of the first diagram. It features the word "ASAP" in a bold, italicized, sans-serif font, with a registered trademark symbol (®) to its upper right. The letters are white with a black outline, set against a red background that represents the vessel wall.

ORDERING INFORMATION

Catalog Number	Description	Box Quantity
ASAP100	ASAP Aspiration Catheter & Kit	1
ASAPLP	ASAPLP Catheter & Kit	1
MMSF200	MiniStop Waste Basin with Lid and Filter	20

U.S. Patent Nos. 6,077,698; 6,278,018; 6,603,040; 6,924,390; 7,138,541



Merit Medical Systems, Inc. • 1600 West Merit Parkway • South Jordan, Utah 84095 • 1-801-253-1600 • 1-800-35-MERIT

Merit Medical Europe, Middle East, and Africa (EMEA) • Amerikalaan 42, 6199 AE Maastricht-Airport • The Netherlands • Tel: +31 43 358 82 22

Merit Medical Ireland Ltd. • Parkmore Business Park West • Galway, Ireland • +353 (0) 91 703 733

Free phone for specific country: Austria 0800 295 374 • Belgium 0800 72 906 (Dutch) 0800 73 172 (French) • Denmark 80 88 00 24 • France 0800 91 60 30

Finland 0800 770 586 • Germany 0800 182 0871 • Ireland (Republic) 1800 553 163 • Italy 800 897 005 • Luxembourg 8002 25 22 • Netherlands 0800 022 81 84

Norway 800 11629 • Sweden 020 792 445 • UK 0800 973 115

www.merit.com

401951010/A ID 053013