

# EU DECLARATION OF CONFORMITY

## MANUFACTURER

Rexxam Co., Ltd. Kagawa Factory  
958 Ikeuchi, Konan-cho,  
Takamatsu-shi, Kagawa-ken,  
761-1494 Japan

SRN : JP-MF-000004361

## AUTHORIZED EUROPEAN REPRESENTATIVE

Rexxam Co., Ltd. Netherlands Office  
Steenhovenweg 5, 5708HN Helmond, The Netherlands

SRN : NL-AR-000002643

## MEDICAL DEVICE

**Model Name** : SLIT LAMP RS-500  
**Classification** : Class I (based on the rule 13 of MDR Annex VIII)  
**Serial number** : effective from 70DS0015 (Main Unit)  
: effective from 70DZ0015 (Power Supply Box)  
**Intended Purpose** : A transparent body such as the cornea or crystalline lens of the eyeball is irradiated with slit light from an oblique direction to perform light cutting, Microscopic abnormalities are raised by the diffused light from the tissue, and magnification observation is performed with a stereoscopic microscope.

THE BASIC UDI-DI 45400080000837

This declaration is only valid for the manufactured products with the document released by us of each product.

The undersigned hereby declares that the medical device as specified above conforms with the General Safety and Performance Requirements listed in Annex I of Regulation (EU) 2017/745.

This declaration of conformity is based on the Regulation (EU) 2017/745 Annex IV.

The device that is covered by the present declaration is in conformity with this Regulation and, with any other relevant union legislation that provides for the issuing of an EU declaration of conformity.  
We declare that we are exclusively responsible for this declaration.

Kagawa Factory, 12. May.2022

Place and date of issue



Name, signature and  
position of manufacture  
Development Dept.1 Manager  
Shunji Takashima

RE-A06-B12-B

# RS-5000/500/300

## LED Slit Lamp RS Series



## About Us

### Rex + Max = Rexxam

Rexxam, which means 'the king of the kings', is a respected and reliable brand.

Rexxam is a Japanese company with a celebrated 60 year history. With over 3,000 employees worldwide, Rexxam manufacture a wide range of products for various industries; from factory automation, automobiles and air conditioning systems, to beer and ski boots.

Since 1986, Rexxam has manufactured various high quality products for leading brands in the eye care industry, including SHIN-NIPPON. Rexxam had developed and manufactured products for SHIN-NIPPON since 1993 and in 2014 the company took over the SHIN-NIPPON brand.

We will be bringing high quality ophthalmic equipment to a global market. By combining precision engineering with industry leading innovation and experience in mass production, Rexxam produce unique products to support eye care specialists across the world.

### Quality in vision care, we are Rexxam.



1960  
Foundation of Rexxam

1986  
Rexxam started the development and manufacturing of ophthalmic devices as an OEM supplier

1993  
Rexxam became the main OEM partner for SHIN-NIPPON  
**SHIN-NIPPON**

2014  
Rexxam acquired the SHIN-NIPPON brand  
**SHIN-NIPPON by Rexxam**

2018  
The manufacturer brand Rexxam was inaugurated  
**Rexxam**

# Rexxam

Quality in vision care

Proudly  Made in Japan



RS-5000



RS-500



RS-300



# Rexxam Slit Lamp RS Series

The new Rexxam slit lamp RS series have high quality optics with versatile functionalities and are easy to use. The series consists of three models; RS-5000, RS-500 and RS-300.

The RS-5000 is the high illumination tower model, while the RS-500 is the low illumination compact model with 5 step magnifications. The RS-300 is the primary model with three step magnifications yet with similar level of high performance and features.

The optics of the RS-5000/RS-500/RS-300 are bright and wide with an optimum convergence angle\*(12 degrees). Drastic improvements were applied to the device in order to achieve comfortable and smooth operation by newly designed joystick and ergonomic base design that ensures solid movement to the various desired positions. The newly designed main and background illumination control dials are designed to be easily and intuitively used with one hand. From the wide range of LEDs, we have chosen the LED which is closest to the conventional halogen light colour temperature.

\*convergence angle is for RS-5000/500 only

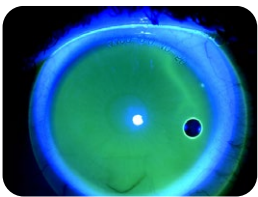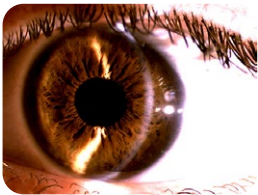

Comfortable to use and high quality optics





The RS-5000 and RS-500 share the same outstanding and high quality features required by slit lamp users. The excellent multicoated optics system allows superb observations and accurate diagnostics. Both the RS-5000 and the RS-500 integrate seamlessly with Rexam's dedicated digital camera DCS-700 (optional accessory) for outstanding digital imaging possibilities with a press of the joystick button.

### LED light source

- Similar colour temperature as Halogen
- Consistent colour temperature at any illumination level
- Low temperature of LED prevents tear film from vaporising
- Cost and energy saving (Long life span and lower energy level)



### Ocular convergence

Comfortable 12 degrees convergence optics for easy stereoscopic fusion



### Magnification

Smooth rotatable drum with 5 step magnification, 6.4x, 10x, 16x, 25x, 40x




### Wide field of view

With wide field optics, we have created an increased field of view. With the well positioned eyepieces, observations are performed in a more natural, comfortable manner. Field of view: 5.9, 9.5, 14.6, 23 and 37mm



### Integrated Yellow filter

The yellow filter is built in as standard, which in combination with the cobalt-blue filter, assures an excellent high contrast examination of the anterior surface.

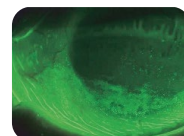

### Slit width and slit length

Rexxam slit lamps offers one of the biggest slit dimensions on the market.

Width : 0-15mm continuous

Length \*: 0.5-15mm continuous

\* RS-5000 only



### Flexible background illumination & Diffuser

The RS-5000/500 are both equipped with background illumination\* and a diffuser that can be easily deployed.

Beside the Cobalt Blue filter, both slit lamps also offer some options as follows: a Heat absorption filter, 50%ND filter\*\* and a Green (Red free) filter

\* background illumination for RS-500 is option

\*\* 50%ND filter not available on RS-500



RS-5000



RS-500

### New combined illumination dials

The illumination dial is a unique feature of the RS-5000/500 combining the regulation of main and background illumination in the same position which allows one hand operation.



### Joystick with shutter button

The newly designed joystick assures smooth movement and incorporates the shutter button for the digital camera DCS-700. The rotation of the joystick supports the weight of the DCS-700 smoothly in vertical movement.



### Ergonomically designed slit lamp base

The base of the RS-5000/500 is newly designed to assure ergonomic grip and smooth movement. The base offers a wide movement range up to 100mm to support the use of fundus lenses. In combination with the use of the DCS-700 digital camera, the base supports automatic Right / Left detection for the camera software user interface.



### New colour

The RS series supports the new colour combination which matches the other Rexxam products, giving your practice an impressive, professional and harmonised look.



# Practical, Economical, LED Slit Lamp RS-300

The RS-300 is the most economical and practical solution in the RS series.

The RS-300 offers 3 step magnification 10x, 16x and 25x in combination with 13.5x eyepieces (13x, 21x and 33x in combination with 18x eyepieces\*). The optical system is of the Galileo type with parallel oculars. The Field of view of 8-20mm is illuminated by an LED light source.

The RS-300 is a high quality slit lamp for primary operations, in any practice.

\*Optional



Applanation Tonometer  
Arm Type  
SL-TM



Applanation Tonometer  
Table Type  
SL-TM



Background Illumination \*\*



Beam Splitter \*



Assistant Scope  
(Monocular)



Assistant Scope  
(Binocular)



C-Mount



— CCD Camera

Measuring Eyepiece\*  
18x Eyepieces \*\*\*



Digital Camera System\*  
DCS-700



\* only for RS-5000/500  
\*\* only for RS500  
\*\*\* only for RS-300





|                       |                                | RS-5000                                                                              | RS-500                                               | RS-300                                                                                                              |
|-----------------------|--------------------------------|--------------------------------------------------------------------------------------|------------------------------------------------------|---------------------------------------------------------------------------------------------------------------------|
| Microscope            | Type                           | Galileo type                                                                         |                                                      |                                                                                                                     |
|                       | Conversion                     | 5 step drum rotation                                                                 |                                                      | 3 step drum rotation                                                                                                |
|                       | Eyepieces                      | 12 degrees                                                                           |                                                      | 0 degree (Parallel)                                                                                                 |
|                       | Magnifications (Field of View) | 12.5x                                                                                |                                                      | 13.5x (18x optional)                                                                                                |
|                       |                                | 6.4x ( 37 mm)<br>10 x ( 23 mm)<br>16 x (14.6 mm)<br>25 x ( 9.5 mm)<br>40 x ( 5.9 mm) |                                                      | 10x (20mm)<br>16x (13mm)<br>25x ( 8mm)<br><br>with 18x optional Eyepieces<br>13x (14mm)<br>21x (10mm)<br>33x ( 6mm) |
|                       | Diopter Adjustment Range       | ±7D                                                                                  |                                                      | ±5D                                                                                                                 |
|                       | Pupillary Distance             | 52 ~ 80mm                                                                            |                                                      | 45 ~ 81mm                                                                                                           |
|                       | Filter                         | Yellow                                                                               |                                                      | —                                                                                                                   |
| Illumination System   | Light Source                   | LED                                                                                  |                                                      |                                                                                                                     |
|                       | Slit Width                     | 0 ~ 15 mm continuous                                                                 |                                                      |                                                                                                                     |
|                       | Slit Length                    | 0.5 ~ 15 mm continuous                                                               | —                                                    |                                                                                                                     |
|                       | Slit Aperture Deameter         | Φ0.2, Φ1, Φ3, Φ5, Φ8, Φ15mm                                                          |                                                      |                                                                                                                     |
|                       | Slit Tilt                      | 5, 10, 15, 20°                                                                       | —                                                    |                                                                                                                     |
|                       | Filters                        | Selectable                                                                           | Cobalt Blue, Heat absorption, Green(Red free), 50%ND |                                                                                                                     |
|                       |                                | Regular Use                                                                          | UV cut                                               |                                                                                                                     |
|                       | Background Illumination        | Yes                                                                                  | optional                                             | —                                                                                                                   |
| Movement Range        | Up - Down                      | 30 mm                                                                                |                                                      |                                                                                                                     |
|                       | Right - Left                   | 100 mm                                                                               |                                                      |                                                                                                                     |
|                       | Forward - Backward             | 100 mm                                                                               |                                                      |                                                                                                                     |
|                       | Fine Movement                  | Φ15 mm (by Joystick)                                                                 |                                                      |                                                                                                                     |
|                       | Chin-rest Up - Down            | 80 mm                                                                                |                                                      |                                                                                                                     |
| Power                 | Power Voltage                  | AC 100 ~ 240 V , 50/60Hz                                                             |                                                      |                                                                                                                     |
|                       | Power Consumption              | 70 VA                                                                                | 50 VA                                                |                                                                                                                     |
| Size                  | Weight                         | 16.5 kg                                                                              | 15.0 kg                                              |                                                                                                                     |
|                       | Dimensions                     | 320mm(W)<br>450mm(D)<br>770mm(H)                                                     | 320mm(W)<br>420mm(D)<br>590mm(H)                     |                                                                                                                     |
| DCS-700 Compatibility |                                | Yes                                                                                  |                                                      | No                                                                                                                  |

## Standard Accessories

- Dust Cover
- Chin-rest Paper
- Test Bar
- Diffuser\*
- Base Plate

\*RS-5000/500 only

Design and specifications are subject to change without prior notice.

### Manufacturer



Rexxam Co.,Ltd.  
Kagawa factory  
958 Ikeuchi, Konan-cho,  
Takamatsu-shi, Kagawa-ken,  
761-1494 Japan

### Contact

Rexxam Co.,Ltd.  
Eye-care Instruments Sales Dept. Tokyo Office  
2-4-2 Kandatsukasa-machi, Chiyoda-ku, Tokyo,  
101-0048 Japan  
TEL:+81-3-6262-9471 FAX:+81-3-6262-9472  
E-mail:eye@rexexam.co.jp  
Website:https://www.rexxam.co.jp

### Distributed by

Proudly Made in Japan

