

Suction Bronchial/Oral

Closed Suction Systems




OptiFlo 72 h safeclean, 56 cm

Closed Endotracheal Suction System safeclean with lockable irrigation chamber, with the possibility of cleaning the catheter without lavage, hence without respiration interruption. The catheter is isolated from the trachea with triple-swivel and vac control.

Application area:	Adult
Total length:	56 cm
Duration of time:	72 h
Material:	Main components: ABS, PP, NBR-white, K-Resin, PU, PVC, PC
Sterile:	Yes

REF	Outer diameter	PU
43.004.57.100	10 CH	60 pcs
43.004.57.120	12 CH	60 pcs
43.004.57.140	14 CH	60 pcs
43.004.57.160	16 CH	60 pcs

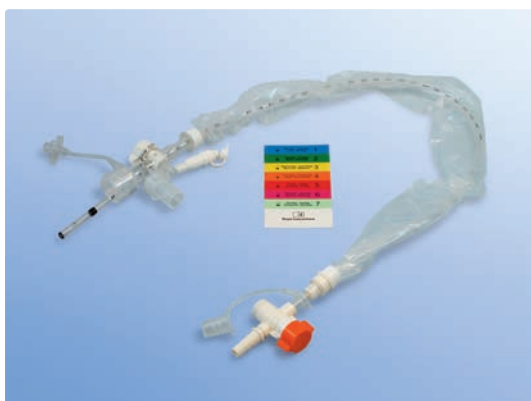

OptiFlo 72 h safeclean, 62 cm

Closed Endotracheal Suction System safeclean with lockable irrigation chamber, with the possibility of cleaning the catheter without lavage, hence without respiration interruption. The catheter is isolated from the trachea with triple-swivel and vac control.

Application area:	Adult
Total length:	62 cm
Duration of time:	72 h
Material:	Main components: ABS, PP, NBR-white, K-Resin, PU, PVC, PC
Sterile:	Yes

REF	Outer diameter	PU
43.004.57.142	14 CH	60 pcs
43.004.57.162	16 CH	60 pcs

Suction Bronchial/Oral

Closed Suction Systems



OptiFlo 72 h safeclean, tracheotomised

Closed Endotracheal Suction System safeclean with lockable irrigation chamber, with the possibility of cleaning the catheter without lavage, hence without respiration interruption. The catheter is isolated from the trachea with triple-swivel and vac control.

Application area:	Adult
Total length:	36 cm
Duration of time:	72 h
Material:	Main components: ABS, PP, NBR-white, K-Resin, PU, PVC, PC
Sterile:	Yes

REF	Outer diameter	PU
43.004.57.121	12 CH	60 pcs
43.004.57.141	14 CH	60 pcs
43.004.57.161	16 CH	60 pcs



OptiFlo Set 48 h 54 cm, catheter mount double swivel, vial

The set contains one OptiFlo closed suction system, 12 vials with 15 ml and catheter mount 15 cm with double swivel, 22F-22M/15F.

Application area:	Adult
Total length (Catheter):	54 cm
Duration of time:	48 h

REF	Outer diameter	PU
43.004.60.143	14 CH	20 pcs
43.004.60.163	16 CH	20 pcs

OptiFlo 72 h safeclean, 56 cm
**REF:**

43.004.57.100, 43.004.57.120, 43.004.57.140, 43.004.57.160

Product group:

Catheters and Probes, sterile

Intended use:

The closed suction system serves as an endotracheal suction catheter to remove airway secretions and/or foreign matter from the airways.

Product description:

The closed safeclean extraction system features a lockable irrigation chamber. This prevents the irrigation fluid from spilling back into the patient when the catheter is closed during irrigation. Furthermore, the system features a straight atraumatic catheter tip, a metered dose aerosol adapter, an irrigation port for cleaning the suction catheter, a PEEP seal, a super-soft catheter sheath, a viewing window and a suction control valve with closure option.

Application area: Adult**Duration of time:** 72 h**Total length:** 56 cm**Sterile:** Yes

REF	Outer diameter	Selling unit quantity	Selling unit GTIN	Packing unit quantity	Packing unit GTIN
43.004.57.100	10 Ch	10 pcs	04034342316457	60 pcs	04034342316464
43.004.57.120	12 Ch	10 pcs	04034342316495	60 pcs	04034342316501
43.004.57.140	14 Ch	10 pcs	04034342316570	60 pcs	04034342316587
43.004.57.160	16 Ch	10 pcs	04034342316693	60 pcs	04034342316709

Materials

Product	ABS, stainless steel, HDPE, K-resin, NBR, PC, PP, PU, PVC, silicone, TPR
Packaging	Tyvek®, cardboard, corrugated cardboard
Latex	No
DEHP/BBP/DBP	No

Additional information

Regulatory basis	Directive 93/42/EEC
Notified Body	TÜV SÜD Product Service GmbH (0123)
Classification	Ila
UMDNS	10-749
eClass	34241302
Manufacturer	P. J. Dahlhausen & Co. GmbH
Single-use Product	Yes
Sterilisation	EO gas
Shelf life	5 Years
Instructions for Use	Yes

OptiFlo 72 h safeclean, 36 cm, tracheotomised
REF: 43.004.57.121, 43.004.57.141, 43.004.57.161

Product group: Catheters and Probes, sterile
Intended use:

The closed suction system serves as an endotracheal suction catheter to remove airway secretions and/or foreign matter from the airways.

Product description:

The closed safeclean extraction system features a lockable irrigation chamber. This prevents the irrigation fluid from spilling back into the patient when the catheter is closed during irrigation. Furthermore, the system features a straight atraumatic catheter tip, a metered dose aerosol adapter, an irrigation port for cleaning the suction catheter, a PEEP seal, a super-soft catheter sheath, a viewing window and a suction control valve with closure option.

Application area: Adult

Duration of time: 72 h

Total length: 36 cm

Sterile: Yes

REF	Outer diameter	Selling unit quantity	Selling unit GTIN	Packing unit quantity	Packing unit GTIN
43.004.57.121	12 Ch	10 pcs	04034342316532	60 pcs	04034342316549
43.004.57.141	14 Ch	10 pcs	04034342316617	60 pcs	04034342316624
43.004.57.161	16 Ch	10 pcs	04034342316730	60 pcs	04034342316747

Materials

Product	ABS, stainless steel, HDPE, K-resin, NBR, PC, PP, PU, PVC, silicone, TPR
Packaging	Tyvek®, cardboard, corrugated cardboard
Latex	No
DEHP/BBP/DBP	No

Additional information

Regulatory basis	Directive 93/42/EEC
Notified Body	TÜV SÜD Product Service GmbH (0123)
Classification	Ila
UMDNS	10-749
eClass	34241302
Manufacturer	P. J. Dahlhausen & Co. GmbH
Single-use Product	Yes
Sterilisation	EO gas
Shelf life	5 Years
Instructions for Use	Yes