



奥美医疗
ALLMED MEDICAL

GAUZE BANDAGE

TECHNICAL PROPERTIES

Content: 100% bleached cotton yarn;

Weft-Warp String number: 40Sx40S/30x30, $\geq 32\text{g/m}^2$;

Width: 10cm, 14cm, other sizes available;

Length: 5m, 7m, Other length can be cut as required;


Shelf life: 5 years

For and on behalf of
ALLMED MEDICAL PRODUCTS CO., LIMITED

Non-Sterile Hydrophilic Gauze Bandage

Susan Liao
Authorized Signature(s)

REF Code	SIZE	BOX in (Roll)	Box Size (cm)
AM-GB10500	10cm x 5m	540	45 x 36.5 x 38 cm
AM-GB14700	14cm x 7m	360	51 x 48.5 x 31 cm

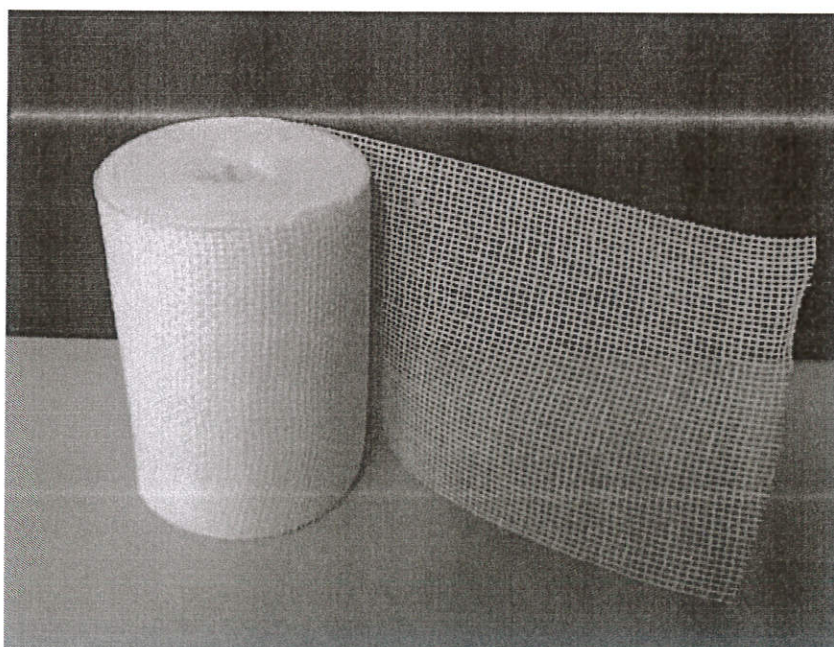

Document No.:	Revision No.:	Date of Revision:	Approved by:
WP-011100-AL-03	R-002	2022.10.12	Susan Liao



奥美医疗
ALLMED MEDICAL

Intended use:

This device is intended to be used for fixation of dressings or tampon on the wound, with protection and barrier function (mechanical) of the injured skin. It provides mild compression and support. MD Class I according to Regulation (EU) 2017/745 of the European Parliament and Council.

Product standards: In compliance with EN 14079.

Certificates: CE, ISO13485.

Labeling: in accordance with Chapter III (23) of Annex I of Regulation 2017/745 of the European Commission.

Cautions:

- Read instructions before use.
- Before usage check the packing completeness, the LOT number and expiry date.
- Use only if package is undamaged and validity not expired;
- For single use.
- Destroy after use.

Properties:

- Non-toxic;
- Hypoallergenic;
- Non-inflammable;
- PHT and DEHP free;
- Latex free;
- Soft and comfortable;
- Low linting;
- High absorbance capacity;
- Easy to cut;
- Free from loose wires;
- Non-Sterile unit packages.

For and on behalf of
ALLMED MEDICAL PRODUCTS CO., LIMITED

Susan Liao
Authorized Signature(s)

Document No.:	Revision No.:	Date of Revision:	Approved by:
WP-011100-AL-03	R-002	2022.10.12	Susan Liao



奥美医疗
ALLMED MEDICAL

Instructions for use:

1. Select the appropriate Gauze Bandage size for your procedure.
2. Examine individual package for integrity of packaging prior to use. Do not use if packaging is broken.
3. Wash hands. (Wear gloves where necessary).
4. Before usage, open the unit package.
5. Apply the bandage in accordance with your facility's protocol.
6. Support the injured area while bandaging (if necessary).
7. Immediately dispose of after single use.

Packing:

- each bandage is rolled and individually packed
- then is packed in strong export carton

Storage:

Protect from sunlight, extreme humidity and heat sources. Not store with organic or chemical solvents. Not store with corrosive gases. Store in dry, cool, good ventilated and clean room, where temperatures between 10°C and 40°C and relative humidity between 30% and 60% are recommended



ALLMED MEDICAL PRODUCTS CO.,LTD.

Address: No.18 Qi Xing Road, Majiadian Town, 443200 Zhijiang City, Hubei Province, P.R.C

For and on behalf of
ALLMED MEDICAL PRODUCTS CO., LIMITED

Susan Liao
.....
Authorized Signature(s)

Document No.:	Revision No.:	Date of Revision:	Approved by:
WP-011100-AL-03	R-002	2022.10.12	Susan Liao