


# BRUCELLOSIS



Brucellosis Antigen  
for the Complement Fixation Test

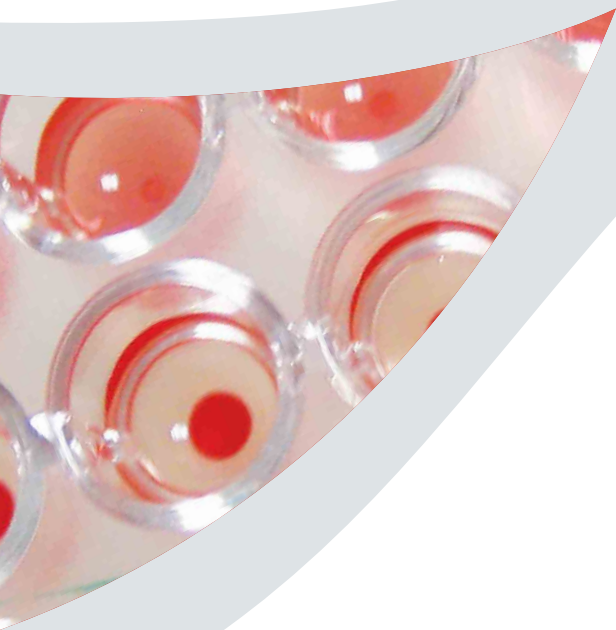

*Brucella abortus* antigen  
for the Complement  
Fixation Test

**ID.vet**  
Innovative Diagnostics

With you at every step

[www.id-vet.com](http://www.id-vet.com)

## Brucellosis Antigen for the Complement Fixation Test

*Brucella abortus* antigen for the Complement Fixation Test.

- Suspension of *Brucella abortus* biovar 1 Weybridge strain No 99, inactivated by heat and phenol.
- Calibrated to give 50% haemolysis with the second international standard serum (OIESS) diluted 1/200

## Ordering information

Product code	AG-BRU
Kit format	100 ml
Reactions	~ 40000 reactions

## References

(1) OIE Terrestrial Manual 2009; 2.4.3; Bovine brucellosis.

## Associated Products

### Brucellosis positive freeze-dried serum

Product Code	MRI-BRU
Format	1 ml vial (freeze-dried)
Description	Freeze-dried highly positive bovine serum to be used as internal reference material for quality control.

### Brucellosis negative freeze-dried serum

Product Code	MRINEG-BRU
Format	1 ml vial (freeze-dried)
Description	Freeze-dried negative bovine serum to be used as internal reference material for quality control.

### Guinea pig complement

Product Code	CPLT-2x5ML
Format	2 x 5 ml vials
Description	Freeze-dried guinea pig complement to be used in any complement fixation test (CFT).

### Veronal buffer

Product Code	VB-4L
Format	4 vials (powder form). Reconstitute each vial with 1L of distilled water
Description	The Veronal Calcium Magnesium buffer is used to dilute sera and reagents for the complement fixation test.