

Sofia[®] PLUS

Soft torqueable catheter *Optimized For Intracranial Access*

Sofia[®]

Distal Access Catheter

Soft torqueable catheter *Optimized For Intracranial Access*

EXTREMELY TRACKABLE

INTRACRANIAL

LARGE LUMEN
SUPPORT CATHETERS

6F


Large inner diameter and enhanced softness for superior distal navigation and procedural versatility




One unit per box, includes shaping mandrel and introducer sheath

PRODUCT CODE	CATHETER SIZE	PROXIMAL OD	DISTAL OD	ID	WORKING LENGTH	DISTAL LENGTH	PROXIMAL LENGTH	DISTAL TIP SHAPE
DA6115ST	6F	0.0825" / 2.1mm	0.0815" / 2.1mm	0.070"	115cm	19cm	96cm	Straight
DA6125ST	6F	0.0825" / 2.1mm	0.0815" / 2.1mm	0.070"	125cm	19cm	106cm	Straight
DA6131ST	6F	0.0825" / 2.1mm	0.0815" / 2.1mm	0.070"	131cm	19cm	112cm	Straight

SOFIA® 6F 115cm

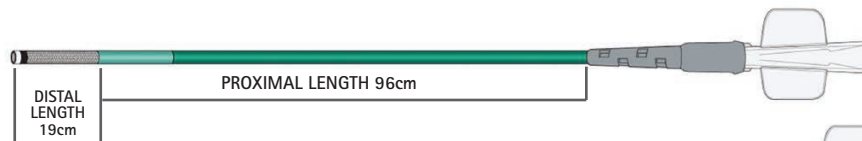

SOFIA® 6F 125cm

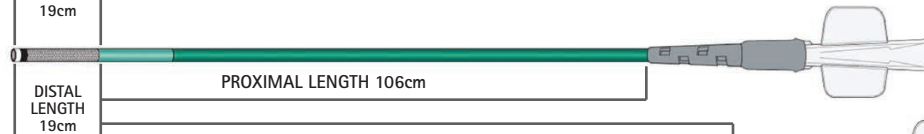

SOFIA® 6F 131cm



Design Features

FEATURES	BENEFITS
0.070" Lumen	Wide Inner Lumen for Capture of Larger Clots*
Exceptionally Soft Distal Tip	Allows Smooth Bypass of the Ophthalmic Artery
Hybrid Braid and Coil Design	Provides Superior 1:1 Push Response
Steam Shapeable Tip and Torqueable Shaft	Ability to Steer Distal Tip Past Difficult Bifurcations
Enhanced Kink Resistance	Maintains Distal and Proximal Lumen Integrity

INDICATIONS FOR USE: The SOFIA® Catheter is indicated for general intravascular use, including the neuro and peripheral vasculature. The SOFIA® Catheter can be used to facilitate introduction of diagnostic or therapeutic agents. The SOFIA® Catheter is not intended for use in coronary arteries.

*with FDA cleared Stentrievors

MicroVention, Inc., Worldwide Headquarters, 1311 Valencia Avenue, Tustin, CA 92780 USA
Phone: 714.247.8000 Customer Service: 800.990.8368



MICROVENTION and SOFIA are registered trademarks of MicroVention, Inc. Refer to Instructions for Use, contraindications and warnings for additional information. © 2016 MicroVention, Inc. MM518 7/16

microvention.com