
# UNIVERSAL SUTURES "HEALING BEYOND COMFORT"



Sutures Mesh Bone Wax

## Needle Shape

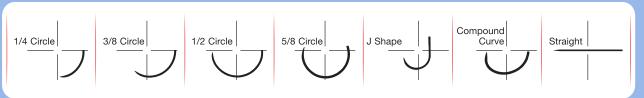

## Needle Type:












International Coding Pattern										
Code 1	Product	Code 2	USP Size	Code 3	Code 4	Needle Curvature	Code 5	Needle Type	Code 6	
1	Braided & Coated Polyglycolic Acid Undyed	Α	3		С	3/8 Circle	Z	Round Body / Taper Point		
2	Braided & Coated Polyglycolic Acid	В	2	H in cm	D	1/2 Circle	Υ	Cutting	in m m	
3	Monofilament Polyamide (Nylon)	С	1		Α	Straight	Х	Reverse Cutting		
4	Catgut Chromic	D	0		E	5/8 Circle	К	Taper Cutting		
5	Black Braided Silk	E	2-0				Т	Trocar Point	LENGTH	
6	Braided Polyester	F	3-0	NGT			R	Blunt Point	8	
7	Catgut Plain	G	4-0	<b>H</b>				Spatula	쁘	
8	Monofilament Polypropylene	Н	5-0	JRE					빌	
9	Monofilament Polydioxanone	ı	6-0	UTUR					NEEDLE	
10	Braided & Coated Polyglactin 910	J	7-0	รเ					۳	
11	Monofilament Polyglecaprone 25	К	8-0							
12	Braided & Coated Polyglactin 910 Undyed									

### Example:

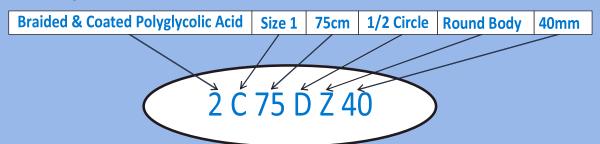

# **UNISIL**

Black Braided Silk Non Absorbable Surgical Suture U.S.P.



**UNISIL – Black Braided Silk** Suture is a non absorbable, sterile, surgical suture composed of an organic protein called fibroin. This protein is derived from the domesticated species Bombyx mori (B. Mori) of the family Bombycidae.

UNISIL sutures are processed to remove the natural waxes and gums. UNISIL suture is dyed black and coated with a special wax mixture. UNISIL Black Braided Silk is also available in its natural color.

UNISIL Virgin silk is available in which the sericin gum is not removed and serves to hold the filaments together. In-house needle manufacturing provides you the better needles with optimum sharpness which finally delivers smooth riding on tissues during suturing.

#### **Indications:**

Black Braided Silk Suture is indicated for use in general soft tissue approximation and/or ligation, including use in cardiovascular, ophthalmic and neurological procedures.

- Highest purified protein from Species B. Mori results no tissue reaction
- · Beeswax coated ultra smooth surface decrease suturing timing followed by no tissue trauma
- Unmatched Knot strength provides satisfaction to user
- USP approved dye increases visibility followed by no reaction

#### **Best Used By:**

- · Orthopaedic Surgeons,
- · General Surgeons,
- · Gynaecologists,
- · Cardiac Surgeons,
- Endoscopic Surgeons

#### **Ultimate Suture For:**

- Skin Closure
- Hard Tissues approximation
- · Tendon Repair
- Muscles approximation

#### Additional Information:

Suture Characteristics : Non Absorbable Surgical Suture

Type : Braided

Material : Braided Silk

Color : Black

Absorption : Non Absorbable

USP Range : 5/0-2

Sterilization : EO (Ethylene Oxide)

Shelf Life : 5 Years

Product Characteristics : Excellent tensile strength, Excellent knot tying, Minimum tissue reaction.

Packaging : In food grade folder packed in a medical grade pouch, with or without needle

Needle Type : Round bodied, cutting edge, taper point, straight, blunt point

MFG. LIC No : MFG/MD/2020/000109

OEM. LIC No : N. Code: KA/DEVICE/MFG/MD/2020/000109

# **UNISIL**

Black Braided Silk Non Absorbable Surgical Suture U.S.P.




Code No.	SUT	URE	NEEDLE DESCRIPTION			
	USP SIZE	LENGTH				
5E75CY26	2-0	75 cm	3/8 Circle Cutting			
5E75CX24	2-0	75 cm	3/8 Circle Reverse Cutting			
5E12X60	2-0	12x60	Without Needle			
5D75CY40	0	75 cm	3/8 Circle Cutting			
5C75CY40	1	75 cm	3/8 Circle Cutting			
5E13X60	2-0	13x60	Without Needle			
5G75CZ16	4-0	75 cm	3/8 Circle Round Bodied 16mm			
5F75DZ30	3-0	75 cm	1/2 Circle Round Bodied 30mm			
5F60	3-0	60 cm	Without Needle			
5D60	0	60 cm	Without Needle			
5F75CX19	3-0	75 cm	3/8 Circle Reverse Cutting 19mm			
5F75CX25	3-0	75 cm	3/8 Circle Reverse Cutting 25mm			
5F75CX30	3-0	75 cm	3/8 Circle Reverse Cutting 30mm			
5F75DX25	3-0	75 cm	1/2 Circle Reverse Cutting 25mm			
5F75DZ26	3-0	75 cm	1/2 Circle Taper Point 26mm			
5F75DX30	3-0	75 cm	1/2 Circle Reverse Cutting 30mm			
5E75DZ22	2-0	75 cm	1/2 Circle Taper Point 22mm			
5E75DZ30	2-0	75 cm	1/2 Circle Taper Point 30mm			
5E75CZ30	2-0	75 cm	3/8 Circle Reverse Cutting 40mm			
5D75DZ26	0	75 cm	1/2 Circle Taper Point 26mm			
5D75DZ35)	0	75 cm	1/2 Circle Taper Point 35mm			
5F45DZ26	3-0	45 cm	1/2 Circle Taper Point 26mm			
5E75DZ25	2-0	75 cm	1/2 Circle Taper Point 25mm			
5D75DZ25	0	75 cm	1/2 Circle Taper Point 25mm			
5F75	3-0	75 cm	Without Needle			
5D60	0	60 cm	Without Needle			









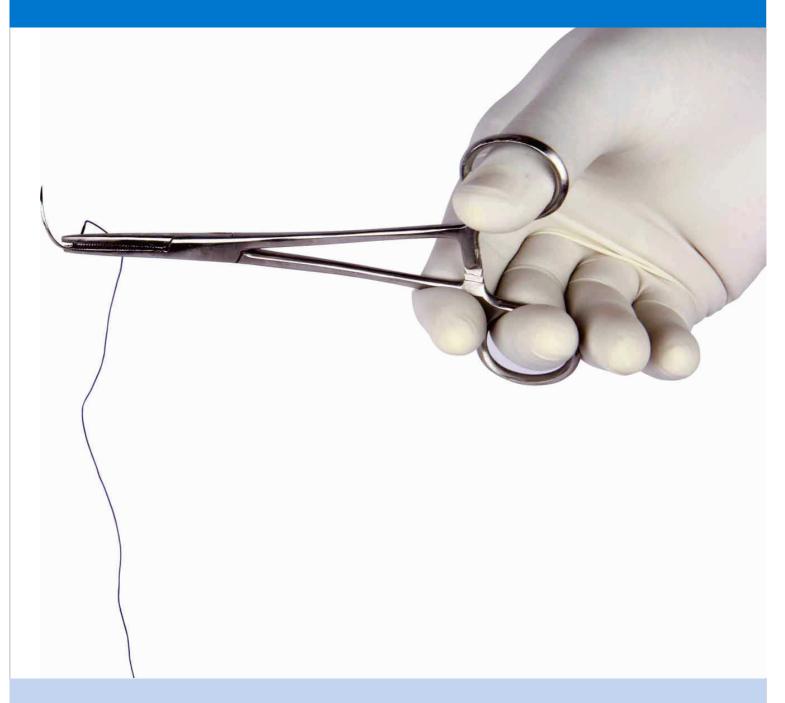

#### REGISTERED OFFICE:

#### UNISUR LIFECARE PVT. LTD.

No. 15/1,2,3 Andhrahalli Main Road, Acharya Indl. Complex, Vishwaneedam Post, Near Anupama School,

BANGALORE - 560091 Karnataka, INDIA

Tel.: +91 9108 990 400 |

 $\hbox{E-Mail: info@universal sutures.com} \mid \hbox{Web.: www.universal sutures.com}$ 

#### CORPORATE OFFICE:

#### UNISUR LIFECARE PVT. LTD.

Unit No. 303, 3rd Floor, Brigade Rubix,

Plot No.: MYS357 in Peenya Plantation, HMT Factory Main Road,

BANGALORE - 560013 Karnataka, INDIA

Tel.: +91 9108 990 400 |

E-Mail: care@universalsutures.com | Web.: www.universalsutures.com



