













# SOFTWARE

Vieworks pursues an "all-in-house design" that offers both hardware and software solutions.

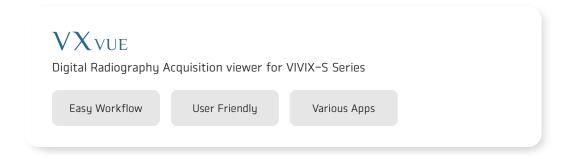

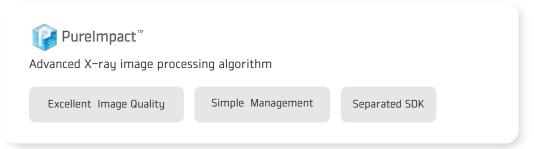

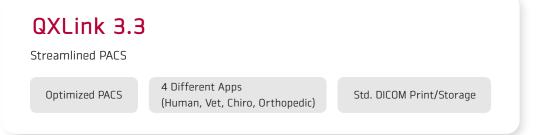

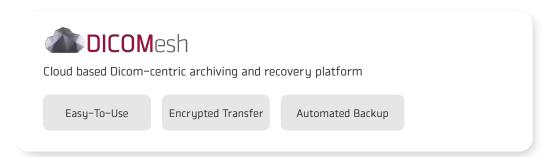



## **VXvue**

### VXvue with PureImpact™

Digital Radiography Acquisition Viewer for VIVIX-S Series

VXvue is an image acquisition program cooperating with VIVIX-S series. New post-processing algorithm, PureImpact™, provides clear images for efficient diagnosis. With various functions and tools, VXvue maximizes work efficiency in different medical environments.



#### Main Features

#### Easy Management

- Main integration systems: Generator, U-arm, collimator, DAP, etc.
- DICOM 3.0 : MWL, Send, Print, MPPS, etc.
- Simple registration with preset auto retrieve
- Exporting the images to external USB/Hard disc

### Advanced Image Processing -PureImpact™

- · New post-processing algorithm for an efficient diagnosis
- Fine details without visual artifact, soft tissue delineation, non-grid chest processing, stable imaging processing within standard radiation dose variation, grid line removal with clear resolution, fine tuning(bright region suppression), wide latitude
- Providing 3 different mode(Soft, Normal, Hard) for individual preference
- Optimized pre-set parameters for each target examination (human, vet, and equine)

## Convenient Usability

- Automated functions including auto-cropping, auto-labeling, and auto-sending
- Touch screen interface
- Customizable UI(font, color, style)
- Statistic function with dose and exam management











Measurement (Free rotation)

### Installation System Requirements

Monitor Resolution (w x h) CPU

Memory

Operating System

1280 X 768

Intel Core i5 2600

Microsoft Windows 8 / 8.1 Pro (64 bit) Microsoft Windows 10 Pro (64 bit)

## QXLink 3.3

#### Streamlined PACS

QXLink 3.3 is a picture archiving system that stores all digital information from diagnostic imaging apparatus. Various measurement tools assist convenient diagnosis.

Users can customize its layout easily with multi-monitor configuration.



#### Main Features

Web-based Access

- Easy management: convenient and unlimited viewer installation with concurrent license
- Easy configuration: convenient configuration with Built-in web server
- Fast delivery: Background and JIT(Just in Time) loading


Multi-monitor Configuration


- Hanging protocol: View with protocol (1–1, 1–1–1, 2–1, etc.) up to 4 monitors and custom view is also available.
- Various screen resolution: VGA to UHD without custom scaling

Advanced Tools with Guide

- More than 60 measurements tools including Chiropractic, Orthopedic, Veterinary
- Useful guide for each tools make it easy to use, minimizing misdiagnosis
- Printing with drag & drop in custom layout is available.
- Full size page for stitched image is available.



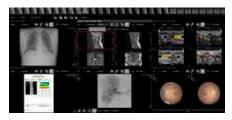

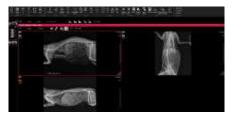

# Installation System Requirements

Monitor Resolution (w x h) CPU Memory Operating System 640 x 480 Intel Core i5 2600 2GB

Microsoft Windows 8 / 8.1 Pro (32 bit / 64 bit)

Microsoft Windows 10 Pro (64 bit)

# DICOMesh Cloud Back-up System

### Synchronized Cloud based archiving and recovery system

Our cloud based Dicom-centric archiving platform provides the ideal solution for your HIPAA & HITECH compliant back-up and recovery needs.

DICOMesh gives the user the freedom in knowing that their back-up system not only provides automated back-up functionality but that the their sensitive data is protected using the latest in HIPAA and HITECH compliant data encryption technology.




#### Main Features

- · Syncronize your local backup to the cloud
- · Able to Query/Retrieve any past uploaded study if local copy is lost
- HIPAA & HITECH compliant solution
- Scalable Storage through structured tiers: 250GB, 500GB, 1TB
- · Perfect DICOM cloud archive
- Finished studies are saved locally to any provided HDD/NAS/Mapped Drive
- Local backups save DICOM & JPEG versions on study completion

# Full Intergration with Vivix Workstations

- Able to Send X-rays (CR, DR, DX) to DICOMesh from VXvue Acquisition Software
- Query/Retrieve lost studies from DICOMesh using QXLink
- · Studies are sent to DICOMesh automatically on study close
- · Annotated and measured in QXLINK? Just send the updated study to DICOMesh

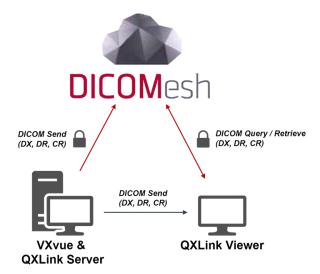

# Installation System Requirements

CPU 1.3 GHz processor minimum

Memory 4GB of free RAM minimum

Storage 3GB free hard drive available space on OS drive

Browser Internet Explorer 11 or higher

Internet Connection TCP/IP Network



a vieworks company

1 Chapin Rd, Unit 3, Pine Brook, NJ 07845 800–456–7800 sales@medlinkimaging.com www.medlinkimaging.com