

# Direct Ophthalmoscopes



# Keeler Direct Ophthalmoscopes

A combination of optical perfection, superb ergonomics and versatile features make Keeler direct ophthalmoscopes the ideal choice. Whatever your requirement from basic examination to detailed diagnosis, there is a Keeler ophthalmoscope for you.

The next couple of pages will guide your choice but first read on to discover the benefits shared by all Keeler ophthalmoscopes as standard.

# Key features

Keeler optical system

The more you see the better the diagnosis. Keeler's wide angle light beam and individual mirror alignment ensures that whatever the patient's pupil size, you will see the clearest view of the retina.

#### Comfort and fit

With curved surfaces and subtle patient side marking, the contoured, curved ophthalmoscope head guarantees an optimal viewing angle and the smooth line surfaces ensure a patient friendly non-threatening instrument.

#### Slimline handles

Keeler hand helds are supplied on a slimline, lightweight robust handle.

#### • Light, tough and built to last

Keeler ophthalmoscopes are built from resilient, lightweight materials to give you many years of trouble free service. High impact resistant covers combined with a strong, rigid chassis protect the carefully aligned optics.

### A clearer view

Controllable light for bright, clear images

The bright white LED, Halogen or Xenon illumination is fully adjustable from zero. This ensures the comfort of your patients and the best possible diagnosis for you.

Reflex free images

Individually controlled optics reduce the effect of corneal reflex whilst minimising shadows and providing clearer retinal visualisation.





#### Apertures, graticules and filters

WIDE ANGLE - Illuminates the largest area of fundus for the best general diagnosis through a dilated pupil.

INTERMEDIATE - Easier access through undilated pupils and in peripheral examination. Ideal for paediatric examination.

MACULAR - For examination of the macular region of the fundus where a larger beam would create excessive pupillary reaction or patient discomfort.

SLIT - To determine retinal elevations and depressions, can also be used to assess anterior chamber depth.

GLAUCOMA - Projects a graticule onto the retina to assess the optic cup/disc ratio as an aid to glaucoma diagnosis.

**FIXATION CROSS** - Projects a graticule onto the retina for assessment of the degree and direction of eccentric fixation.

RED FREE AND FIXATION CROSS - For vessel examination with a fixation cross for assessment of eccentric fixation.

**RED FREE** - Used to examine blood vessels in fine detail. Blood vessels are silhouetted black against a dark green background.

COBALT BLUE - Use with fluorescein dye for the examination of corneal scars and abrasions.

# LED Light Source

Keeler Direct Ophthalmoscopes are now available with a LED light source, giving bright, white illumination and 3 x longer lasting battery life.







3 x Longer lasting battery life

Battery Life	Battery will last 3 x longer with the LED light source than the existing filament bulb at maximum brightness before requiring a re-charge.
Patch Quality	Homogenous, sharp, clear patch.
Colour Temperature	3000K which equates to a neutral coloured patch.
CRI	High-CRI therefore colours are rendered more accurately.
Electronics	Uses the latest in LED technology.
Compatibility	The LED light source is compatible with the Specialist, Professional, Practitioner and Standard Keeler Ophthalmoscope heads and the 3.6V lithium ion rechargeable handle and the 2.8V primary cell handle.



# Specialist

We recognise that illumination is of fundamental importance to good diagnosis. The re-engineered optics give 30% more illumination and because it is longer lasting, that means fewer bulb changes giving greater peace of mind.

- The World's finest Direct Ophthalmoscope
- Lightweight and well balanced

Almost 50% lighter, the Specialist combined with the slimline handle is perfectly balanced.

#### Expanded dioptre range

The most comprehensive lens range available in a direct ophthalmoscope, +44D to -45D lens range guarantees optimal diagnosis in any situation.



# Key features

#### • 30% more illumination

Re-engineered optics give 30% more illumination to get through, for example, cloudy media and cataracts.

#### Slimmer optics

Slimmer and enclosed viewing optics enable you to get closer to the eye providing a wider field of view.

#### • Premium hand made optics

The Morton Lens Track system, made up of individually ground glass lenses ensures an optimal view.

#### New Cobalt Blue filter

The new blue filter can be used in conjunction with fluorescein dye for the detection and examination of corneal scars and abrasions.

#### Pupillometer gauge

The pupillometer disc allows quick and easy estimation of pupil size.

### Part Numbers

#### **Specialist Ophthalmoscopes**

**1132-P-1157** 3.6V Specialist Ophth (rechargeable) in case

1132-P-1149 2.8V Specialist Ophth in case

**1132-P-1181** 3.6V LED Specialist Ophth (rechargable) in case

1132-P-1173 2.8V LED Specialist Ophth in case

#### Accessories

1011-P-7034 Spec/Vista halogen bulb 3.6V Pkt2

1011-P-7042 Spec/Vista halogen bulb 2.8V Pkt2

1011-P-7229 Spec/Vista LED bulb 2.8V/3.6V Pk1



# Professional

The Professional Ophthalmoscope has been designed to ensure ultimate control. Lenses, graticules and filters can all be introduced at the touch of a fingertip without the need to pull away from the patient.

Positive action maywheel and comprehensive lens range

The positive action, easy access, maywheel not only allows you to dial up the required lens quickly, but also gives you total control of where you are in the lens range making each exam quicker and easier.

The extensive range of +29D to -30D in single dioptre steps ensures easy examination of all ocular structures and covers a wide range of patient/user prescriptions. A quick step swing of +/- 20D helps speed up diagnosis.







### Key features

Choice of illumination

The Professional now has a choice of Xenon or LED illumination for a brighter and whiter long lasting light. Giving you 3 x longer battery life.

Swing over filter

The red free filter for enhanced vessel examination can be used in conjunction with any of the 6 diagnostic beams.

Brow rest

The integral soft brow rest allows for easy positioning and protects your lenses if you are a spectacle wearer.



### Part Numbers

#### **Professional Ophthalmoscopes**

1128-P-1014 3.6V Professional Ophth in case
1128-P-1015 3.6V LED Professional Ophth in case
1128-P-1013 2.8V Professional Ophth in case
1128-P-1008 2.8V LED Professional Ophth in case

#### Accessories

 1011-P-7114
 Std/Pract/Prof Xenon Bulb 3.6V Pkt2

 1011-P-7106
 Std/Pract/Prof Xenon Bulb 2.8V Pkt2

 1011-P-7210
 Std/Pract/Prof LED Bulb 3.6/2.8V Pkt1



# Professional Retinoscope

Optimised for outstanding performance with either a Streak or Spot halogen beam with a simple bulb change. Plus the LED Professional Streak Retinoscope is now available!

• Neutralisation check (Patented)

The only Retinoscope allowing rapid confirmation of neutralisation. You can change the beam from divergent to convergent enabling you to cross check the accuracy of the neutral point, guaranteeing the assessment is correct every time.

Precision magnetic control (Patented)

Improved 360° access to a single control for both rotation and convergence adjustments. The precision movements are achieved by use of rare earth magnet technology.



### Key features

Choice of illumination

Choose between 600 lux of brilliant halogen illumination or the LED Retinoscope. Whichever you choose, the Professional guarantees you enough light to refract even the most difficult eyes.

LED Professional Streak Retinoscope
 Gives you better contrast, image quality and
 longer battery life.

Two position aperture

A simple finger movement allows a choice of 4mm or 1.7mm aperture to optimise brightness or resolution of the retinal reflex.

### Part Numbers

#### Professional Retinoscopes

**1302-P-1014** 3.6V Professional Streak LED Ret and case **1302-P-1011** 3.6V Professional Streak Retinoscope in case

1302-P-1009 3.6V Professional Spot Retinoscope in case 1302-P-1013 2.8V Professional Streak LED Ret and case

1302-P-1010 2.8V Professional Streak LED Ret and case

1302-P-1008 2.8v Professional Spot Retinoscope in case

#### Accesso

1305-P-7000 Streak Retinoscope LED blister pack

#### **Professional Combi Retinoscope Bulbs**

1013-P-7009 Vista/Professional Streak bulb 3.6V Pkt2

1013-P-7011 Professional Spot bulb 3.6V Pkt2

1013-P-7010 Professional Spot bulb 2.8V Pkt2

#### **Vista Retinoscope Bulbs**

1013-P-7009 Vista/Professional Streak bulb 3.6V Pkt2

1013-P-7007 Vista Spot halogen bulb 3.6V Pkt2

1013-P-7006 Vista Spot halogen bulb 2.8V Pkt2

 $^{\circ}$ 

# Ophthalmoscope/Retinoscope Sets

You can purchase your Ophthalmoscope and Retinoscope in a slimline double handle case and add in a Lithium Double Charger.



Specialist Ophthalmoscope/Retinoscope Set



### Part Numbers

#### Specialist Ophthalmoscope/Retinoscope Sets

1427-P-104 3.6V Streak LED Ret and Specialist Ophth
1427-P-1063 Slimline 3.6V Streak Ret and Specialist Ophth
1427-P-1065 Slimline 3.6V Spot Ret and Specialist Ophth
1427-P-1103 2.8V Streak LED Ret and Specialist Ophth
1427-P-1062 Slimline 2.8V Streak Ret and Specialist Ophth
1427-P-1064 Slimline 2.8V Spot Ret and Specialist Ophth
1427-P-1116 3.6V Specialist & Combi LED only Student Set
1427-P-1119 3.6V Specialist Student Set LED only Volk

#### Professional Ophthalmoscope/Retinoscope sets

1427-P-102 3.6V Streak LED Ret and Professional Ophth
1427-P-1059 3.6V Streak/Professional Set with 2 handles
1427-P-1061 Slimline 3.6V Spot Ret and Professional Ophth
1427-P-1101 2.8V Streak LED Ret and Professional Ophth
1427-P-1058 Slimline 2.8V Streak Ret & Professional Ophth
1427-P-1060 Slimline 2.8V Spot Ret & Professional Ophth
1427-P-1071 3.6V Magnetic Streak Ret/Professional Ophth and handle
1427-P-1117 3.6V Professional LED only Student Set

1427-P-1118 3.6V Professional Student Set LED only Volk

**Lithium Double Charger** 

1941-P-1368 Lithium Double Charger

1999-P-7308 Double Charger Wall Mount Kit

# Ophthalmoscope/Retinoscope Wall Units

Keeler's anti-theft designed Ophthalmoscopes and Retinoscopes have a lockable head and handle fitting so you'll always find them where you left them.



If you require a wall unit you will need to order the transformer unit which comes with the 2 x handles and each individual head.

# Key features

Tangle-free

The strong, coiled, tangle-free cable permits instrument use up to 3 metres from the transformer.

Durable

The robust construction is ideal for busy working conditions.

Automated

Automatic illumination when instruments are removed from the cradle. With LED confirmation.



### Part Numbers

#### Wall Units

1964-P-2001 Double Cord Handle Power Supply
1132-P-1050 3.6V Specialist Ophth Head and Bulb
1128-P-1004 3.6V Professional Ophth in Carton
1132-P-1042 2.8V Specialist Ophth Head and Bulb
1128-P-1003 2.8V Professional Ophth in Carton
1302-P-1012 Professional Streak LED
Retinoscope Head
1302-P-1007 3.6V Streak Head in Carton
1302-P-1005 3.6V Spot Magnetic Retinoscope
1302-P-1004 2.8V Spot Magnetic Retinoscope
1302-P-1004 2.8V Spot Magnetic Retinoscope



# Accessories

#### Colour Coded Grips

Add a bit of colour to your diagnostic area with colour coded grips. Available in four colours these new grips are not only fun but functional too!

Available individually or as an assorted pack of four, choose from pink, green, blue or black.





Near Point Cards for Retinoscopy

	is compare			
yesterday	step turn	arrows	direction	prien
accurate	glasses test	ing shor	tcut s	ummaris
child er	TOF		effec	t ever
remind	list		matu	red gi
price syst	em		flag	flows
confirm				note loft
justification	television	machine	optical	observ
recognise	examine	environ	ment	calcular
179-1209				E3Keck





#### Keeler Slimline Handles

Spare handles for your instruments available in 3.6V Lithium rechargeable (red base) or 2.8V Battery (blue base) options.





### Part Numbers

#### **Colour Coded Grips**

1901-P-7028	Pink
1901-P-7036	Green
1901-P-7044	Blue
EP29-05365	Black
1901-P-7052	Pack of four assorted
	Colour Coded Grips

#### **Near Point Cards for Retinoscopy**

1321-P-7021 Near Point Cards for Retinoscopy

#### **Keeler Slimline Handles**

**1911-P-1201** 3.6V Slimline Handle **1901-P-1101** 2.8V Slimline Handle

# Chargers and Upgrade Kits

A two hour fast charge and over 2 hours continuous use means your instrument is always ready for use. Choose between the Lithium-ion Double or Mini Chargers.

#### • Lithium-ion Double Charger

The compact charger takes two handles and fits neatly into your immediate testing environment, allowing easy access to your hand held instruments. It can be desk or wall mounted.

#### • Lithium-ion Mini Charger

Can be connected quickly and simply to the base of your Keeler handle, making it the perfect choice for domiciliary visits or practitioners on the move. Your instruments can still be used whilst the battery is charging.



## Key features

#### • Over 2 hours of continuous use

Ideal for long examinations especially while first learning with your instrument.

#### • LED displays

LEDs indicate when battery life is low, charging or fully charged.

#### Upgrade Kits available

Keep your existing heads and upgrade to Lithium-ion.

The Double Lithium Upgrade Kit (pictured below) includes 2 x handles and a Lithium-ion Double Charger.

The Mini Charger Upgrade Kit includes 1 x handle and a Lithium-ion Mini Charger.



### Part Numbers

#### **Lithium-ion Chargers**

1941-P-1368 Lithium Double Charger1941-P-1341 Lithium Mini Charger

#### **Upgrade Kits**

1911-P-1228 Lithium Duo Charger and 2 Slimline Handles

1911-P-1236 Slimline Handle and Mini Charger

#### **Lithium Battery**

EP39-18918 3.6V Lithium Battery

# What's in the box?

#### Specialist

• Specialist Ophthalmoscope 1132-P-1157 with Halogen Bulb 3.6V SPECIALIST • Lithium Battery OPHTHALMOSCOPE • Spare Halogen Bulb (RECHARGEABLE) IN CASE • User Instructions • Specialist Ophthalmoscope 1132-P-1149 with Halogen Bulb 2.8V SPECIALIST • Spare Halogen Bulb OPHTHALMOSCOPE IN CASE • User Instructions 1132-P-1181 • Specialist Ophthalmoscope with LED Bulb 3.6V LED SPECIALIST Lithium Battery OPHTHALMOSCOPE (RECHARGABLE) Case IN CASE • User Instructions • Specialist Ophthalmoscope 1132-P-1173 2.8V LED SPECIALIST with LED Bulb **OPHTHALMOSCOPE** • Case IN CASE User Instructions

#### Professional

1128-P-1014 3.6V PROFESSIONAL OPHTHALMOSCOPE IN CASE	<ul> <li>Professional Ophthalmoscope with Halogen Bulb</li> <li>Lithium Battery</li> <li>Spare Halogen Bulb</li> <li>Case</li> <li>User Instructions</li> </ul>
1128-P-1015 3.6V LED PROFESSIONAL OPHTHALMOSCOPE IN CASE	<ul> <li>Professional Ophthalmoscope with LED Bulb</li> <li>Lithium Battery</li> <li>Case</li> <li>User Instructions</li> </ul>
1128-P-1013 2.8V PROFESSIONAL OPHTHALMOSCOPE IN CASE	<ul> <li>Professional Ophthalmoscope with Halogen Bulb</li> <li>Spare Halogen Bulb</li> <li>Case</li> <li>User Instructions</li> </ul>
1128-P-1008 2.8V LED PROFESSIONAL OPHTH IN CASE	Professional Ophthalmoscope with LED Bulb     Case     User Instructions

#### Retinoscope

1302-P-1014 3.6V PROFESSIONAL STREAK LED RETINOSCOPE AND CASE	<ul> <li>LED Streak Retinoscope</li> <li>Lithium Battery</li> <li>Eyebrow Extension</li> <li>Case</li> <li>User Instructions</li> </ul>
1302-P-1011 3.6V PROFESSIONAL STREAK RETINOSCOPE IN CASE	<ul> <li>Halogen Streak Retinoscope</li> <li>Spare Halogen Bulb</li> <li>Lithium Battery</li> <li>Eyebrow Extension</li> <li>Case</li> <li>User Instructions</li> </ul>
1302-P-1009 3.6V PROFESSIONAL SPOT RETINOSCOPE IN CASE	<ul> <li>Halogen Spot Retinoscope</li> <li>Spare Halogen Bulb</li> <li>Lithium Battery</li> <li>Eyebrow Extension</li> <li>Case</li> <li>User Instructions</li> </ul>
1302-P-1013 2.8V PROFESSIONAL STREAK LED RETINOSCOPE AND CASE	LED Streak Retinoscope     Eyebrow Extension     Case     User Instructions
1302-P-1010 2.8V PROFESSIONAL STREAK RETINOSCOPE IN CASE	<ul><li> Halogen Streak Retinoscope</li><li> Spare Halogen Bulb</li><li> Eyebrow Extension</li><li> Case</li><li> User Instructions</li></ul>
1302-P-1008 2.8V PROFESSIONAL SPOT RETINOSCOPE IN CASE	<ul> <li>Halogen Spot Retinoscope</li> <li>Spare Halogen Bulb</li> <li>Eyebrow Extension</li> <li>Case</li> <li>User Instructions</li> </ul>

### Specialist Ophthalmoscope/ Retinoscope Sets

1427-P-1104 3.6V STREAK LED RETINOSCOPE AND SPECIALIST OPHTHALMOSCOPE	<ul> <li>LED Streak Retinoscope</li> <li>LED Specialist Ophthalmoscope</li> <li>Lithium Battery (x2)</li> <li>Eyebrow Extension</li> <li>Double Handheld Case</li> <li>User Instructions</li> </ul>
1427-P-1063 SLIMLINE 3.6V STREAK RETINOSCOPE AND SPECIALIST OPHTHALMOSCOPE	<ul> <li>Halogen Streak Retinoscope</li> <li>Halogen Specialist</li> <li>Ophthalmoscope</li> <li>Spare Halogen Streak Bulb</li> <li>Spare Halogen Specialist Bulb</li> <li>Lithium Battery (x2)</li> <li>Eyebrow Extension</li> <li>Double Handheld Case</li> <li>User Instructions</li> </ul>
1427-P-1065 SLIMLINE 3.6V SPOT RETINOSCOPE AND SPECIALIST OPHTHALMOSCOPE	<ul> <li>Halogen Spot Retinoscope</li> <li>Halogen Specialist</li> <li>Ophthalmoscope</li> <li>Spare Halogen Spot Bulb</li> <li>Spare Halogen Specialist Bulb</li> <li>Lithium Battery (x2)</li> <li>Eyebrow Extension</li> <li>Double Handheld Case</li> <li>User Instructions</li> </ul>
1427-P-1103 2.8V STREAK LED RETINOSCOPE AND SPECIALIST OPHTHALMOSCOPE	LED Streak Retinoscope     LED Specialist Ophthalmoscope     Eyebrow Extension     Double Handheld Case     User Instructions
1427-P-1062 SLIMLINE 2.8V STREAK RETINOSCOPE AND SPECIALIST OPHTHALMOSCOPE	<ul> <li>Halogen Streak Retinoscope</li> <li>Halogen Specialist</li> <li>Ophthalmoscope</li> <li>Spare Halogen Streak Bulb</li> <li>Spare Halogen Specialist Bulb</li> <li>Eyebrow Extension</li> <li>Double Handheld Case</li> <li>User Instructions</li> </ul>
1427-P-1064 SLIMLINE 2.8V SPOT RETINOSCOPE AND SPECIALIST OPHTHALMOSCOPE	<ul> <li>Halogen Spot Retinoscope</li> <li>Halogen Specialist</li> <li>Ophthalmoscope</li> <li>Spare Halogen Spot Bulb</li> <li>Spare Halogen Specialist Bulb</li> <li>Eyebrow Extension</li> <li>Double Handheld Case</li> <li>User Instructions</li> </ul>

### Student LED Specialist Ophthalmoscope/ Retinoscope Sets

	• LED Streak Retinoscope
	• LED Specialist Ophthalmoscope
	Spare Halogen Spot Bulb
	• Lithium Battery (x2)
	Eyebrow Extension
	User Instructions
	Dynamic Targets
1427-P-1116	Lithium Double Charger
3.6V SPECIALIST AND	Lithium Mini Charger
COMBI LED ONLY	Keeler Lens Cloth
STUDENT SET	Student How to Use
	Ophthalmoscope and
	Retinoscope Book
	Romanes Occulder Rule
	<ul> <li>Fixation Bar Budgie Stick</li> </ul>
	Keeler Reusable Pen Torch
	Double Handheld Soft Case
	Student City Bag
	LED Streak Retinoscope
	·
	·
	• LED Specialist Ophthalmoscope
	LED Specialist Ophthalmoscope     Spare Halogen Spot Bulb
	<ul><li>LED Specialist Ophthalmoscope</li><li>Spare Halogen Spot Bulb</li><li>Volk Digital Wild Field Lens</li></ul>
	<ul> <li>LED Specialist Ophthalmoscope</li> <li>Spare Halogen Spot Bulb</li> <li>Volk Digital Wild Field Lens</li> <li>Lithium Battery (x2)</li> </ul>
	<ul> <li>LED Specialist Ophthalmoscope</li> <li>Spare Halogen Spot Bulb</li> <li>Volk Digital Wild Field Lens</li> <li>Lithium Battery (x2)</li> <li>Eyebrow Extension</li> </ul>
1427-P-1119	<ul> <li>LED Specialist Ophthalmoscope</li> <li>Spare Halogen Spot Bulb</li> <li>Volk Digital Wild Field Lens</li> <li>Lithium Battery (x2)</li> <li>Eyebrow Extension</li> <li>User Instructions</li> </ul>
3.6V SPECIALIST	<ul> <li>LED Specialist Ophthalmoscope</li> <li>Spare Halogen Spot Bulb</li> <li>Volk Digital Wild Field Lens</li> <li>Lithium Battery (x2)</li> <li>Eyebrow Extension</li> <li>User Instructions</li> <li>Dynamic Targets</li> </ul>
3.6V SPECIALIST STUDENT SET LED	<ul> <li>LED Specialist Ophthalmoscope</li> <li>Spare Halogen Spot Bulb</li> <li>Volk Digital Wild Field Lens</li> <li>Lithium Battery (x2)</li> <li>Eyebrow Extension</li> <li>User Instructions</li> <li>Dynamic Targets</li> <li>Lithium Double Charger</li> </ul>
3.6V SPECIALIST	<ul> <li>LED Specialist Ophthalmoscope</li> <li>Spare Halogen Spot Bulb</li> <li>Volk Digital Wild Field Lens</li> <li>Lithium Battery (x2)</li> <li>Eyebrow Extension</li> <li>User Instructions</li> <li>Dynamic Targets</li> <li>Lithium Double Charger</li> <li>Lithium Mini Charger</li> </ul>
3.6V SPECIALIST STUDENT SET LED	LED Specialist Ophthalmoscope     Spare Halogen Spot Bulb     Volk Digital Wild Field Lens     Lithium Battery (x2)     Eyebrow Extension     User Instructions     Dynamic Targets     Lithium Double Charger     Lithium Mini Charger     Keeler Lens Cloth
3.6V SPECIALIST STUDENT SET LED	<ul> <li>LED Specialist Ophthalmoscope</li> <li>Spare Halogen Spot Bulb</li> <li>Volk Digital Wild Field Lens</li> <li>Lithium Battery (x2)</li> <li>Eyebrow Extension</li> <li>User Instructions</li> <li>Dynamic Targets</li> <li>Lithium Double Charger</li> <li>Lithium Mini Charger</li> <li>Keeler Lens Cloth</li> <li>Student How to use</li> </ul>
3.6V SPECIALIST STUDENT SET LED	<ul> <li>LED Specialist Ophthalmoscope</li> <li>Spare Halogen Spot Bulb</li> <li>Volk Digital Wild Field Lens</li> <li>Lithium Battery (x2)</li> <li>Eyebrow Extension</li> <li>User Instructions</li> <li>Dynamic Targets</li> <li>Lithium Double Charger</li> <li>Lithium Mini Charger</li> <li>Keeler Lens Cloth</li> <li>Student How to use</li> <li>Ophthalmoscope and</li> </ul>
3.6V SPECIALIST STUDENT SET LED	LED Specialist Ophthalmoscope Spare Halogen Spot Bulb Volk Digital Wild Field Lens Lithium Battery (x2) Eyebrow Extension User Instructions Dynamic Targets Lithium Double Charger Lithium Mini Charger Keeler Lens Cloth Student How to use Ophthalmoscope and Retinoscope Book Romanes Occulder Rule
3.6V SPECIALIST STUDENT SET LED	LED Specialist Ophthalmoscope     Spare Halogen Spot Bulb     Volk Digital Wild Field Lens     Lithium Battery (x2)     Eyebrow Extension     User Instructions     Dynamic Targets     Lithium Double Charger     Lithium Mini Charger     Keeler Lens Cloth     Student How to use     Ophthalmoscope and     Retinoscope Book
3.6V SPECIALIST STUDENT SET LED	LED Specialist Ophthalmoscope     Spare Halogen Spot Bulb     Volk Digital Wild Field Lens     Lithium Battery (x2)     Eyebrow Extension     User Instructions     Dynamic Targets     Lithium Double Charger     Lithium Mini Charger     Keeler Lens Cloth     Student How to use     Ophthalmoscope and     Retinoscope Book     Romanes Occulder Rule     Fixation Bar Budgie Stick

# What's in the box?

#### Professional Ophthalmoscope/ Retinoscope Sets

• LED Streak Retinoscope 1427-P-1102 LED Professional Ophthalmoscope 3.6V STREAK LED • Lithium Battery (x2) **RETINOSCOPE AND**  Eyebrow Extension **PROFESSIONAL** • Double Handheld Case OPHTHALMOSCOPE User Instructions • Halogen Streak Retinoscope • Halogen Professional Ophthalmoscope 1427-P-1059 • Spare Halogen Streak Bulb 3.6V STREAK/ • Spare Halogen Specialist Bulb PROFESSIONAL SET • Lithium Battery (x2) WITH 2 HANDLES • Eyebrow Extension • Double Handheld Case User Instructions • Halogen Spot Retinoscope • Halogen Professional 1427-P-1061 SLIMLINE Ophthalmoscope 3.6V SPOT • Spare Halogen Spot Bulb RETINOSCOPE AND • Spare Halogen Specialist Bulb **PROFESSIONAL** • Lithium Battery (x2) OPHTHALMOSCOPE Eyebrow Extension • Double Handheld Case User Instructions • LED Streak Retinoscope 1427-P-1101 LED Professional 2.8V STREAK LED Ophthalmoscope **RETINOSCOPE AND**  Eyebrow Extension **PROFESSIONAL** • Double Handheld Case **OPHTHALMOSCOPE**  User Instructions • Halogen Streak Retinoscope • Halogen Professional 1427-P-1058 Ophthalmoscope **SLIMLINE 2.8V** • Spare Halogen Streak Bulb STREAK RETINOSCOPE • Spare Halogen Specialist Bulb AND PROFESSIONAL • Eyebrow Extension OPHTHALMOSCOPE • Double Handheld Case • User Instructions • Halogen Spot Retinoscope • Halogen Professional 1427-P-1060 Ophthalmoscope **SLIMLINE 2.8V SPOT** • Spare Halogen Spot Bulb RETINOSCOPE AND • Spare Halogen Specialist Bulb **PROFESSIONAL** • Evebrow Extension **OPHTHALMOSCOPE** • Double Handheld Case

User Instructions

• Halogen Professional Ophth 1427-P-1071 • 1 Handle shared between 3.6V MAGNETIC Retinoscope and Ophth STREAK RETINOSCOPE/ • Spare Halogen Streak Bulb **PROFESSIONAL** • Spare Halogen Specialist Bulb **OPHTHALMOSCOPE** • Lithium Battery AND HANDLE • Eyebrow Extension • Double Handheld Case User Instructions

• Halogen Streak Retinoscope

#### Student LED Professional Ophthalmoscope/ Retinoscope Sets

• Spare Halogen Spot Bulb • Lithium Battery (x2) • Eyebrow Extension User Instructions • Dynamic Targets • Lithium Double Charger 3.6V PROFESSIONAL · Lithium Mini Charger • Keeler Lens Cloth • Student How to Use Ophth and Retinoscope Book • Romanes Occulder Rule • Fixation Bar Budgie Stick • Keeler Reusable Pen Torch • Double Handheld Soft Case • Student City Bag • LED Streak Retinoscope

• LED Streak Retinoscope

• LED Professional Ophthalmoscope

1427-P-1118 3.6V PROFESSIONAL STUDENT SET LED **ONLY VOLK** 

1427-P-1117

STUDENT SET

**LED ONLY** 

• Dynamic Targets • Lithium Double Charger • Lithium Mini Charger • Keeler Lens Cloth • Student How to Use Ophth and

• LED Professional Ophthalmoscope

• Spare Halogen Spot Bulb

• Lithium Battery (x2)

• Eyebrow Extension

User Instructions

• Volk Digital Wild Field Lens

Retinoscope Book • Romanes Occulder Rule • Fixation Bar Budgie Stick • Keeler Reusable Pen Torch • Double Handheld Soft Case • Student City Bag

#### Wall Units

• Handheld Wall Mount • 2 x Corded Handles 1964-P-2001 • 1.5mm Hexagon Wrench DOUBLE CORD • Ball Ended Hex Wrench 2mm HANDLE • Charger **POWER SUPPLY** • Wall Bracket, Wall Plugs & Screws • User Instructions 1132-P-1050 • Specialist Ophthalmoscope 3.6V SPECIALIST Head with Halogen Bulb OPHTHALMOSCOPE User Instructions **HEAD AND BULB** 1128-P-1004 3.6V PROFESSIONAL • Professional Ophthalmoscope **OPHTHALMOSCOPE** Head with Halogen Bulb IN CARTON 1302-P-1012 • Retinoscope Head with LED **PROFESSIONAL** Streak Bulb STREAK LED • Eyebrow Extension **RETINOSCOPE HEAD** • User Instructions • Retinoscope Head with 1302-P-1007 Streak Bulb 3.6V STREAK HEAD • Eyebrow Extension IN CARTON • User Instructions • Retinoscope Head with 1302-P-1005 3.6V SPOT Spot Bulb MAGNETIC • Eyebrow Extension RETINOSCOPE • User Instructions

#### **Finoff Transilluminator**

2903-P-7000

• Finoff Transilluminator Head 2.8V **FINOFF** TRANSILLUMINATOR • User Instructions 2.8V 2903-P-7001 • Finoff Transilluminator Head 3.6V FINOFF TRANSILLUMINATOR • User Instructions 3.6V

# Specifications

#### **Electrical Ratings**

The ophthalmoscope/retinoscope/otoscope, the power supply (EP29-32777) with its charging dock (1941-P-5289 and 1941-P-5326) together constitute a Medical Electrical System as defined in EN 60601-1.

INPUT MAINS DATA	100-240V – 50/60Hz
POWER SUPPLY RATING	12V : 2.5amps
Operation	
RUN TIME	Maximum 15 minutes ON Minimum 10 minutes OFF

#### Instrument Heads and Handles

#### **INPUT VOLTAGE (DC)** 3V 2xAA Alkaline Batteries - BLUE 3.75V Lithium Ion Rechargeable Battery - RED (EP39-18918) 3.65V NiMH Rechargeable Battery - BLACK (1919-P-7149)

#### Dimensions and Weight

Difficisions and we	EIGHT
SPECIALIST	260 x 35 x 40mm (HxWxD) Weight 270g
PROFESSIONAL	212 x 43 x 28mm (HxWxD) Weight 260g
RETINOSCOPE	240 x 32 x 40mm (HxWxD) Weight 260g
FINOFF TRANSILLUMINATOR	203 x 44 x 29mm (HxWxD) Weight 217g
Classification	
EQUIPMENT	Class II
CLASSIFICATION	Type B protection against shock

### Overview and intended use

The Keeler Ophthalmoscope is indicated for examining the posterior segment of the eye referred to as the fundus, in order to aid in the screening and diagnosis of retinal pathology including, but not limited to, diseases such as cataracts, papilloedema, glaucomatous disc cupping, diabetic retinopathy, hypertensive retinopathy and retinal detachments. When set to a high power and magnification, it can also be used to examine the anterior segment of the eye, which includes the eyelids, cornea, sclera, conjunctiva, iris, aqueous, crystalline lens and anterior vitreous.

The Keeler Retinoscope is indicated for the objective assessment of refractive state of the eye. By observing the retinal red reflex, it is also a means to providing information regarding the visual system, for example media and lens opacities, significant ocular aberrations and accommodative status.

The Keeler Finoff Transilluminator is indicated for the transillumination of tissue, and, in particular, scleral transillumination.



#### More Keeler Instruments

With over 100 years of industry experience, Keeler has become a household name among optometrists and ophthalmologists worldwide.

View our full range at www.keeler.co.uk | www.keelerusa.com

### Got a question? - Call us

UK: +44 (0)1753 857177 India: +91 99 3031 1090

USA: +1 610 353 4350 Brazil: +55 11 4302 6053

China: +86-18512119109









Device Name(s):

# EC Declaration of Conformity- Direct Instruments

Revision: 11
Review Date: 05/02/2020

Manufacturer: Keeler Ltd.

Clewer Hill Rd., Windsor, Berkshire, SL4 4AA

Address: United Kingdom

a. Ophthalmoscope, battery powered

b. Ophthalmoscope, line powered

c. Retinoscope, battery powered

d. Retinoscope, line powered

e. Otoscope, direct

Device Classification: Class I

a. 46786 Type P

b. 46787 Type P

GMDN Classification Code(s) c. 32712 Type P

d. 46784 Type P

e. 12849 Type P

I, the undersigned, hereby declare that the medical devices specified above and bearing CE marking, conform to the applicable provisions of the Medical Devices Directive (MDD) 93/42/EEC and the transpositions of this amended Directive into national law of the countries of sale.

Keeler hereby certifies that the part numbers listed are in compliance with EU Directive 2011/65/EU (RoHS)

This declaration of conformity is issued under the Sole Responsibility of the manufacturer on the basis of:

- Quality Management System Certificate No: LRQ 0943149/E Issued by Lloyds Register Quality Assurance
- ii. Technical File demonstrating conformity with the essential requirements of this Directive and maintained according to MDD Annex VII for devices classified as Class I under MDD Annex IX
- iii. Each device being released to market in accordance with records certifying compliance with applicable product specifications and procedures.

Full Name: Moria Zigman Signature:

Position: Senior Regulatory Affairs Date: 05/02/2019

Officer

Keeler Ltd Clewer Hill Road Windsor Berkshire SL4 4AA Tel: +44 (0) 1753 857177 Email: <a href="mailto:info@keeler.co.uk">info@keeler.co.uk</a>
Fax: +44 1753 827145 (Customer Services) Fax: +44 (0) 1753 830247 (Manufacturing) Website: <a href="mailto:www.keeler.co.uk">www.keeler.co.uk</a>

Registered in England No.408759 Registered office: Clewer Hill Road Windsor Berkshire SL4 4AA VAT No. GB 349 0761 40

QM05 Reference 1 DofC Template Issue 3

bu Zegnon



Revision: 11

Review Date: 05/02/2020

#### Ophthalmoscope

Part Number	<b>Device Description</b>
1102-P-1092	JAZZ POCKET OPHTHALMOSCOPE 2.8V
1102-P-1105	JAZZ LED POCKET OPHTHALMOSCOPE 2.8V
1102-P-1113	JAZZ ULTRA POCKET OPHTHALMOSCOPE 2.8V
1126-P-1005	2.8V STANDARD OPHTH IN CASE
1126-P-1006	3.6V STANDARD OPHTH IN CASE
1126-P-1011	2.8V STANDARD LED OPHTH IN CASE
1126-P-1012	3.6V STANDARD LED OPHTH IN CASE
1126-P-5005	STANDARD OPHTH + LED MODULE
1127-P-1005	2.8V PRACTITIONER OPHTH IN CASE
1127-P-1006	3.6V PRACTITIONER OPHTH IN CASE
1127-P-1008	2.8v PRACTITIONER LED OPHTH IN CASE
1127-P-1009	3.6v PRACTITIONER LED OPHTH IN CASE
1127-P-5003	PRACTITIONER OPHTH + LED MODULE
1128-P-1003	2.8V PROFESSIONAL OPHTH IN CARTON
1128-P-1004	3.6V PROFESSIONAL OPHTH IN CARTON
1128-P-1007	3.6V PROFESSIONAL LED OPHTH IN CARTON
1128-P-1008	2.8V LED PROFESSIONAL OPHTHALMOSCOPE
1128-P-1013	2.8V PROFESSIONAL OPHTH IN CASE
1128-P-1014	3.6V PROFESSIONAL OPHTH IN CASE
1128-P-1014-500	3.6V Professional Ophth in Case
1128-P-1015	3.6V LED PROFESSIONAL OPHTHALMOSCOPE
1128-P-5008	PROFESSIONAL OPHTH + LED MODULE
1130-P-1000	WIDE ANGLE OPHTH IN CARTON
1131-P-5228	SPECIALIST OPHTH + LED MODULE
1132-P-1149	2.8V SPECIALIST OPHTH IN CASE
1132-P-1157	3.6V SPECIALIST OPHTH (RECHARGEABLE)
1132-P-1157-500	3.6V Specialist Ophth in Case
1132-P-1173	2.8V LED SPECIALIST OPHTHALMOSCOPE
1132-P-1181	3.6V LED SPECIALIST OPHTHALMOSCOPE
1132-P-1202	LED SPECIALIST OPHTH HD ONLY
2903-P-7000	2.8V TRANSILLUMINATOR IN CARTON
2903-P-7001	3.6V TRANSILLUMINATOR IN CARTON

#### Retinoscope

Part Number	Device Description
1302-P-1004	2.8V SPOT MAGNETIC RET
1302-P-1005	3.6V SPOT MAGNETIC RET
1302-P-1006	2.8V STREAK MAGNETIC RET
1302-P-1008	2.8V PROFESSIONAL SPOT RETINOSCOPE IN CASE
1302-P-1009	3.6V PROFESSIONAL SPOT RETINOSCOPE IN CASE
1302-P-1010	2.8V PROFESSIONAL STREAK RETINOSCOPE IN CASE
1302-P-1011	3.6V PROFESSIONAL STREAK RETINOSCOPE IN CASE

QM05 Reference 1 DofC Template



Revision: 11
Review Date: 05/02/2020

1302-P-1013	2.8V PROFESSIONAL STREAK LED RETINOSCOPE & CASE
1302-P-1014	3.6V PROFESSIONAL STREAK LED RETINOSCOPE & CASE
1321-P-7021	DYNAMIC TARGETS FOR RETINOSCOPE
1302-P-1016	3.6V LED STREAK RETINOSCOPE IN CASE
1302-P-1011-500	3.6V Professional Streak Retinoscope in Case
1302-P-1014-500	3.6V Professional Streak LED Retinoscope in Case

#### Otoscope

Part Number	Device Description
1503-P-1034	2.8V POCKET OTO IN CASE
1502-P-1091	2.8V DELUXE OTO IN CASE
1502-P-1104	3.6V DELUXE OTO IN CASE
1513-P-1062	2.8V PRACTITIONER OTO
1514-P-1136	SLIMLINE 2.8V FIBRE-OPTIC OTO IN CASE
1514-P-1144	SLIMLINE 3.6V FIBRE-OPTIC OTO IN CASE
1514-P-1160	JAZZ LED POCKET OTOSCOPE 2.8V
1514-P-1187	JAZZ ULTRA POCKET OTOSCOPE 2.8V
1516-P-1025	2.8V STANDARD OTO IN CARTON
1516-P-5069	KEELER 2.8V STANDARD OTOSCOPE
1516-P-5077	KEELER 3.6V STANDARD OTOSCOPE
1518-P-1001	SLIMLINE 2.8V STANDARD OTO
1518-P-1002	SLIMLINE 3.6V STANDARD OTO
1519-P-1001	SLIMLINE 2.8V PRACTITIONER OTO
1519-P-1002	3.6V PRACTITIONER OTO IN CASE

#### Sets

Part Number	Device Description
1102-P-1041	POCKET SET 2.8V
1427-P-1053	PROFESSIONAL & COMBI RET STUDENT SET
1427-P-1054	SPECIALIST & COMBI RET SET
1427-P-1058	SLIMLINE 2.8V STREAK RET & PROFESSIONAL OPHTH
1427-P-1059	3.6V STREAK/PROFESSIONAL SET WITH 2 HANDLES
1427-P-1060	SLIMLINE 2.8V SPOT RET & PROFESSIONAL OPHTH
1427-P-1061	SLIMLINE 3.6V SPOT RET & PROFESSIONAL OPHTH
1427-P-1062	SLIMLINE 2.8V STREAK RET & SPECIALIST OPHTH
1427-P-1063	SLIMLINE 3.6V STREAK RET & SPECIALIST OPHTH
1427-P-1064	SLIMLINE 2.8V SPOT RET & SPECIALIST OPHTH
1427-P-1065	SLIMLINE 3.6V SPOT RET & SPECIALIST OPHTH
1427-P-1066	3.6V STREAK RET & SPECIALIST OPHTH
1427-P-1071	3.6V MAGNETIC STREAK RET/PROFESSIONAL OPHTH & HANDLE
1427-P-1082	AUSTRALIAN SPECIALIST SET
1427-P-1083	AUSTRALIAN SPECIALIST SET
1427-P-1097	2.8V INDIAN SET - OPHTH +LED RET
1427-P-1098	3.6V STREAK RET & PROFESSIONAL OPHTH WITH 2 HANDLES



Revision: 11
Review Date: 05/02/2020

1427-P-1101	2.8V STREAK LED RETINOSCOPE & PROFESSIONAL OPHTHALMOSCOPE
1427-P-1102	3.6V STREAK LED RETINOSCOPE & PROFESSIONAL
	OPHTHALMOSCOPE
1427-P-1103	2.8V STREAK LED RETINOSCOPE & SPECIALIST
	OPHTHALMOSCOPE
1427-P-1104	3.6V STREAK LED RETINOSCOPE & SPECIALIST
	OPHTHALMOSCOPE
1427-P-1107	3.6V STREAK LED RET & PROFESSIONAL OPHTH & 1 HANDLE
1427-P-1108	AUSTRALIAN LED SPECIALIST SET
1427-P-1109	AUSTRALIAN LED PROFESSIONAL SET
1427-P-1113	2.8V LED STREAK RET & PROF SINGLE
1427-P-1114	2.8V INDIAN OPHTH & LED RET WITH 2 HANDLES
1427-P-1115	3.6V STREAK RET LED PROF OPHTH WITH 2 HANDLES
1427-P-1125	2.8V STK LED RET & PRACT OPHTH WITH 1 HANDLE
1427-P-1126	3.6V STK LED RET & PRACT OPHTH WITH 1 HANDLE
1427-P-1128	2.8V STK LED RET & PROF OPHTH - INDIA
1427-P-1129	3.6V STK LED RET & PROF OPHTH - INDIA
1501-P-7117	MAG LENS/ATC FOR STD+PKT OTO
1514-P-1179	JAZZ BULB OTOSCOPE SET
1702-P-1037	POCKET SET 2.8V
1727-P-1009	VISTA FIBRE-OPTIC SET DRY 2.8V HAL
1728-P-1350	GENMED RECH SET A 230V UK
1728-P-1473	GENMED CORD SET A 230V UK
1729-P-1018	SLIMLINE 2.8V STANDARD SET
1729-P-1033	SLIMLINE 2.8V PROFESSIONAL/FIBRE-OPTIC OTO SET
1,23 1 1033	SERVICINE 2.07 I NOI ESSIONAL, I IDRE OF HE OF OSET

#### **Desk Sets**

Part Number	Device Description
1729-P-1028	3.6V PRACTITIONER DESK SET
1729-P-1029	3.6V STANDARD DESK SET
1729-P-1030	3.6V PRACTITIONER/FIBRE-OPTIC DESK SET
1729-P-1057	3.6V PRACTITIONER LED / FO DESK SET
1729-P-1055	3.6V STANDARD LED DESK SET
1729-P-1056	3.6V PRACTITIONER LED DESK SET

#### **Diagnostic Set**

Part Number	Device Description
1702-P-1088	JAZZ POCKET DIAGNOSTIC SET 2.8V
1702-P-1096	JAZZ LED POCKET DIAGNOSTIC SET 2.8V
1702-P-1109	JAZZ ULTRA POCKET DIAGNOSTIC SET 2.8V
1729-P-1019	SLIMLINE 2.8V PRACTITIONER DIAGNOSTIC SET
1729-P-1020	SLIMLINE 2.8V PRACTITIONER/FIBRE-OPTIC DIAGNOSTIC SET
1729-P-1022	SLIMLINE 3.6V STANDARD DIAGNOSTIC SET



Revision: 11
Review Date: 05/02/2020

1729-P-1023	SLIMLINE 3.6V PRACTITIONER DIAGNOSTIC SET
1729-P-1024	SLIMLINE 3.6V PRACTITIONER/FIBRE-OPTIC DIAGNOSTIC SET
1729-P-1024-500	3.6V Pract / FO Diagnostic Set
1729-P-1041	2.8V STANDARD/FIBRE-OPTIC DIAGNOSTIC SET
1729-P-1044	2.8V STANDARD VET DIAGNOSTIC SET
1729-P-1046	2.8V PROFESSIONAL/PRACTITIONER DIAGNOSTIC SET
1729-P-1047	2.8V PRACTITIONER VET DIAGNOSTIC SET
1729-P-1048	2.8V STANDARD LED OPHTH / OTO DIAGNOSTIC SET
1729-P-1049	3.6V STANDARD LED OPHTH / OTO DIAGNOSTIC SET
1729-P-1050	2.8V PRACTITIONER LED OPHTH / OTO DIAGNOSTIC SET
1729-P-1052	2.8V STANDARD LED/FIBRE-OPTIC DIAGNOSTIC SET
1729-P-1053	2.8V PRACTITIONER LED/FIBRE-OPTIC DIAGNOSTIC SET
1729-P-1054	3.6V PRACTITIONER LED/FIBRE-OPTIC DIAGNOSTIC SET
1964-P-2002	3.6V PRACTITIONER DIAGNOSTIC SET
1964-P-2003	3.6V STANDARD DIAGNOSTIC SET
1964-P-2004	3.6V PRACTITIONER FIBRE-OPTIC DIAGNOSTIC SET

#### **Student Sets**

Part Number	Device Description
1427-P-1057	PROFESSIONAL STUDENT KIT, SLIMLINE HANDLE (USA)
1427-P-1067	3.6V SPECIALIST/STREAK STUDENT SET
1427-P-1068	3.6V PROFESSIONAL STUDENT SET
1427-P-1084	US SPECIALIST STUDENT KIT
1427-P-1085	STUDENT SET (SOUTH AMERICA)
1427-P-1086	2.8V SPOT RET STUDENT SET (ITALY)
1427-P-1087	2.8V STREAK RET STUDENT SET (ITALY)
1427-P-1088	2.8V STREAK OPHTH STUDENT SET (ITALY)
1427-P-1089	3.6V STREAK/PROFESSIONAL STUDENT SET
1427-P-1090	3.6V SPECIALIST/STREAK STUDENT SET
1427-P-1091	3.6V STREAK/PROF STUDENT SET (HOLLAND)
1427-P-1092	2.8V STUDENT SET (INDIA)
1427-P-1093	2.8V STUDENT SET (TAIWAN)
1427-P-1094	3.6V STUDENT SET (TAIWAN)
1427-P-1095	3.6V STUDENT SET CW
1427-P-1100	US PROFESSIONAL LED STUDENT KIT
1427-P-1100-500	US Professional LED Student Set
1427-P-1105	SPECIALIST & COMBI LED STUDENT SET
1427-P-1110	3.6V PROFESSIONAL LED STUDENT SET
1427-P-1116	3.6V SPECIALIST & COMBI LED ONLY STUDENT SET
1427-P-1117	3.6V PROFESSIONAL LED ONLY STUDENT SET
1427-P-1118	3.6V PROFESSIONAL STUDENT SET LED ONLY
1427-P-1119	3.6V SPECIALIST STUDENT SET LED ONLY
1427-P-1124	3.6V PROF/STREAK LED STUDENT HOLLAND
1427-P-1127	US PROFESSIONAL COMPLETE LED STUDENT SET



Revision: 11
Review Date: 05/02/2020

1427-P-1127-500 US Professional Complete LED Student Set

#### **Vet Sets**

Part Number	Device Description
1729-P-1038	3.6V STANDARD VET DIAGNOSTIC SET
1729-P-1039	3.6V PRACTITIONER VET DIAGNOSTIC SET
1729-P-1043	3.6V STANDARD VET DESK SET
1729-P-1045	3.6V PRACTITIONER VET DESK SET
1964-P-2013	3.6V CORDED STANDARD VET SET
1964-P-2014	3.6V CORDED PRACTITIONER VET SET

#### **Accessories**

Part Number	Device Description
1011-P-5522	POCKET OPHTHALMOSCOPE LED ASSEMBLY
1011-P-5530	JAZZ OPHTHALMOSCOPE BULB REF 211.04.01
1011-P-5610	OPHTH LED ASSY
1011-P-5637	SPECIALIST LED ASSY
1011-P-7034	SPEC/VISTA HAL BULB 3.6V PKT2
1011-P-7042	SPEC/VISTA HAL BULB 2.8V PKT2
1011-P-7050	POCKET OPHTH HAL BULB 2.8V PKT2
1011-P-7106	STD/PRACT/PROF XE BULB 2.8V PKT2
1011-P-7114	STD/PRACT/PROF XE BULB 3.6V PKT2
1011-P-7122	JAZZ ULTRA OPHTHALMOSCOPE LED
1011-P-7202	JAZZ OPHTHALMOSCOPE BULB PKT2
1011-P-7210	1x OPHTH LED BLISTER PACK
1011-P-7229	1x SPECIALIST LED BLISTER PACK
1013-P-7000	OLD SPOT RET HAL BULB 2.8V PKT2
1013-P-7001	OLD SPOT RET HAL BULB 3.6V PKT2
1013-P-7002	OLD STREAK RET HAL BULB 2.8V PKT2
1013-P-7003	OLD STREAK RET HAL BULB 3.6V PKT2
1013-P-7004	PURVIS STREAK HAL BULB 2.8V PKT2
1013-P-7005	PURVIS STREAK HAL BULB 3.6V PKT2
1013-P-7006	VISTA SPOT HAL BULB 2.8V PKT2
1013-P-7007	VISTA SPOT HAL BULB 3.6V PKT2
1013-P-7008	VISTA/PROF STREAK BULB 2.8V PKT2
1013-P-7009	VISTA/PROF STREAK BULB 3.6V PKT2
1013-P-7010	PROF SPOT BULB 2.8V PKT2
1013-P-7011	PROF SPOT BULB 3.6V PKT2
1015-P-5298	POCKET OTOSCOPE LED ASSEMBLY
1015-P-7023	STD OTOSCOPE HAL BULB 3.6V PKT2
1015-P-7031	STD OTOSCOPE HAL BULB 2.8V PKT2
1015-P-7058	FO OTOSCOPE HAL BULB 3.6V PKT2
1015-P-7066	FO OTOSCOPE HAL BULB 2.8V PKT2
1015-P-7074	JAZZ ULTRA OTOSCOPE SPARE LED

QM05 Reference 1 DofC Template



Revision: 11
Review Date: 05/02/2020

POKT OPH BALL + SPRING CAN 1101-P-5032 2.8V POCKET OPHTHALMOSCOPE HEAD & BULB 1102-P-1033 KEELER POCKET LED OPHTH HEAD + LED 1102-P-5077 JAZZ ULTRA OPHTH LFD HFAD 1102-P-5085 2.8V STANDARD OPHTH HEAD 1126-P-1002 3.6V STANDARD OPHTH HEAD & BULB 1126-P-1004 STANDARD LED OPHTH HEAD IN CART 1126-P-1013 1126-P-8015 PROFESSIONAL OPHTH MAYWHEEL REPAIR KIT 1126-P-8016 **OPHTH FIXINGS KIT** 2.8V PRACTITIONER OPHTH HEAD IN CARTON 1127-P-1002 3.6V PRACTITIONER OPHTH HEAD & BULB 1127-P-1004 1127-P-1010 PRACTITIONER LED HEAD IN CARTON 1128-P-1017 LED PROFESSIONAL OPHTHALMOSCOPE - HEAD ONLY WIDE ANGLE OPHTH HEAD IN CARTON 1130-P-5004 1132-P-1042 2.8V SPECIALIST OPHTH HEAD & BULB 1132-P-1050 3.6V SPECIALIST OPHTH HEAD & BULB 1132-P-1165 3.6V SPECIALIST OPHTH HEAD & LED 1302-P-1007 3.6V STREAK HEAD IN CARTON PROFESSIONAL STREAK LED RETINOSCOPE HEAD 1302-P-1012 1305-P-7000 STREAK RET LED BLISTER PACK BAG OF 500 DISPOSABLE SPECULUM NO.2 1501-P-5170 BAG OF 500 DISPOSABLE SPECULUM NO.4 1501-P-5189 BAG OF 500 NO.1 (2.5MM) SPECULA 1501-P-5197 BAG OF 500 NO.3 (4.5MM) SPECULA 1501-P-5218 1501-P-5306 BAG OF 100 NO.1 SPECULA (2.5MM) BAG OF 100 NO.2 SPECULA (3.5MM) 1501-P-5314 1501-P-5322 BAG OF 100 NO.3 SPECULA (4.5MM) 1501-P-5330 BAG OF 100 NO.4 SPECULA (5.5MM) 1512-P-7008 SPECULUM NO.1 1512-P-7016 SPECULUM NO.2 1512-P-7024 SPECULUM NO.3 1512-P-7032 SPECULUM NO.4 1512-P-7040 SPECULUM NO.5 VISTA AIR TIGHT CHAMBER-ATC MAC 1512-P-7067 1512-P-7083 AIRTIGHT CHAMBER/PNEUMATIC TESTING 1513-P-1054 3.6V PRACTITIONER OTO HEAD & BULB TUBE OF DISPOSABLE SPECULA 1513-P-5020 1513-P-7026 DISPOSABLE SPECULA (NO ADAPTOR) 1513-P-7034 VIEWING LENS & MNT-VISTA OTO 2.8V VISTA FIBRE-OPTIC OTO HEAD & BULB 1514-P-1064 1514-P-1072 3.6V VISTA FIBRE-OPTIC HEAD & BULB 1514-P-5110 KEELER POCKET FIBRE OPTIC OTOSCOPE HEAD 1514-P-5137 JAZZ ULTRA OTOSCOPE HEAD 1514-P-7001 FIBRE OPTIC/STANDARD INSUFFLATOR ADAPTOR

A HALMA COMPANY



Revision: 11
Review Date: 05/02/2020

PNEUMATIC TESTING ADAPTER 1514-P-7028 1514-P-7036 JAZZ 2MM REUSABLE SPECULUM (PK of 10) JAZZ 2.5MM REUSABLE SPECULUM (PK of 10) 1514-P-7044 JAZZ 3MM REUSABLE SPECULUM (PK of 10) 1514-P-7052 1514-P-7060 JAZZ 4MM REUSABLE SPECULUM (PK of 10) JAZZ 5MM REUSABLE SPECULUM (PK of 10) 1514-P-7079 1514-P-7087 JAZZ 2MM SPECULUM (PK of 100) 1514-P-7095 JAZZ 2.5MM SPECULUM (PK of 100) 1514-P-7108 JAZZ 3MM SPECULUM (PK of 100) 1514-P-7116 JAZZ 4MM SPECULUM (PK of 100) 1514-P-7124 JAZZ 5MM SPECULUM (PK of 100) 1514-P-7132 JAZZ 2.5MM/4.0MM DISPOSABLE SPECULUM (PK - 5 of EACH) 1516-P-1033 3.6V VISTA STANDARD OTO HEAD & BULB 1517-P-1019 2.8V POCKET OTO HEAD 1599-P-7018 **BULB EXTRACTOR** 1599-P-7229 SET OF PERMANENT SPECULA 1599-P-7245 INSUFFLATOR BULB FOR ALL OTO KEELER TORCH - WITH REPLACEBLE BATTERY 1602-P-1064 **KEELER SILVER PENLIGHTS - PACK OF 6** 1602-P-1072 PENLIGHT, SILVER 1602-P-1080 RECHARGE DIAGNOSTIC KIT A-240V 1728-P-1000 **RECHARGE DIAGNOSTIC KIT B-240V** 1728-P-1051 RECHARGE DIAGNOSTIC KIT C-240V 1728-P-1115 **CORD HANDLE DIAG KIT A-240V** 1728-P-1174 CORD HANDLE DIAG KIT B-240V 1728-P-1238 **CORD HANDLE DIAG KIT C-240V** 1728-P-1289 1901-P-1101 2.8V SLIMLINE HANDLE 1901-P-5356 **POCKET HANDLE** POCKET HANDLE DRY CELL BOTTOM CAP 1901-P-5372 2.8V POCKET HANDLE 1901-P-5380 1901-P-5444 BLUE BASE CAP, 2.8V SLIMLINE HANDLE (AA) 1901-P-5452 RED BASE CAP, 3.6V SLIMLINE HANDLE (LITHIUM) 1901-P-5460 **KEELER JAZZ HANDLE** 1901-P-5479 JAZZ SHORT HANDLE ASSY SLIMLINE HANDLE SLEEVE - PINK 1901-P-7028 SLIMLINE HANDLE SLEEVE - GREEN 1901-P-7036 1901-P-7044 SLIMLINE HANDLE SLEEVE - BLUE 1901-P-7052 SLIMLINE HANDLE SLEEVE ASSORTED COLOURS SLIMLINE LITHIUM HANDLE 1911-P-1201 1911-P-1201-500 Slimline Rechargeable Handle + Battery LITHIUM DUO CHARGER & 2 SLIMLINE HANDLES 1911-P-1228 SLIMLINE HANDLE & MINI CHARGER 1911-P-1236 1911-P-1244-500 Rechargeable Handle + Refract Stand Adapt 1919-P-7093 3.6V LITHIUM HANDLE UPGRADE KIT

QM05 Reference 1 DofC Template

Issue 3



Revision: 11

Review Date: 05/02/2020

1941-P-1341	LITHIUM MINI CHARGER
1941-P-1368	LITHIUM DOUBLE CHARGER
1941-P-1384	WIRELESS CHARGER
1941-P-5289	LITHIUM DOUBLE CHARGER ONLY
1941-P-5326	LITHIUM MINI CHARGER ASSY
1964-P-2000	SPECULA DISPENSER UNIT
1964-P-2001	DOUBLE CORD HANDLE POWER SUPPLY
1964-P-2006	GENMED/DISPOSE-A-SPEC & BRACKET
1964-P-2020	INTEGRATED DIAGNOSTIC WALL KIT - ACCOSON
1964-P-2021	INTEGRATED DIAGNOSTIC WALL KIT
1964-P-7001	SLIMLINE CORDED WALL BRACKET KIT
1999-P-5660	SLIMLINE HANDLE SUB-ASSY
1999-P-7308	DOUBLE CHARGER WALL MOUNT KIT
2202-P-1000	Keeler Retinoscopy Rack
EP29-05365	SLIMLINE HANDLE SLEEVE - BLACK
EP29-32777	Power Supply
EP39-18918	3.6V LITHIUM BATTERY

EC Certificate Production Quality Assurance System: Certificate GB20/965237



The management system of

### **Keeler Ltd**

Clewer Hill Road, Windsor, Berkshire, SL4 4AA, UK

has been assessed and certified as meeting the requirements of

### Directive 93/42/EEC

on medical devices, Annex V

For the following products
Applanation Tonometer to aid diagnosis
and measurement of intraocular pressure

Disposable Applanation Tonometer Cone for use in testing of intraocular pressure

Disposable Cryo Probe for use with Ophthalmic Surgery Devices

Where the above scope includes class IIb or class III medical device(s), a valid EC Type Examination

Certificate according to Annex III is a mandatory requirement for each device
in addition to this certificate to place that device on the market.

This certificate is valid from 28 February 2020 until 30 September 2022 and remains valid subject to satisfactory surveillance audits.

Issue 1. Certified since 08 September 2009 and first certified by SGS Belgium NV since 28 February 2020

Certification is based on reports numbered GB/PC/ 240569

Authorised by


SGS Belgium NV, Notified Body 1639

SGS House Noorderlaan 87 2030 Anwerp Belgium t +32 (0)3 545-48-48 f +32 (0)3 545-48-49 www.sgs.com

LPMD5008 - Certificate CE1639 AnnexV EN rev. 01

Page 1 of 1



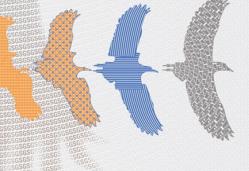

This document is a Web version of SGS certificate for electronic use exclusively. It shall only be available by clicking on SGS Certification Mark which has been posted on Your website. It shall not be printed in anyway. This document is copyright protected. No content or appearance may be reproduced without the express written permission of SGS. Any misuse, alteration, forgery or falsification is unlawful.

Certificate GB19/964765



The management system of

### **Keeler Ltd**

Clewer Hill Road, Windsor, Berkshire, SL4 4AA, UK

has been assessed and certified as meeting the requirements of

ISO 13485:2016 EN ISO 13485:2016

For the following activities

Design, Development and manufacture, distribution & service of Optometric and Ophthalmic medical devices and equipment.

This certificate is valid from 02 December 2019 until 31 October 2022 and remains valid subject to satisfactory surveillance audits.

Re certification audit due before 04 September 2022

Issue 1. Certified since 10 September 1995

Authorised by



SGS United Kingdom Ltd Rossmore Business Park Ellesmere Port Cheshire CH65 3EN UK t +44 (0)151 350-6666 f +44 (0)151 350-6600 www.sgs.com

HC SGS 13485 2016 0118

Page 1 of 1











This document is issued by the Company subject to its General Conditions of Certification Services accessible at www.sgs.com/terms\_and\_conditions.htm. Attention is drawn to the limitations of liability, indemnification and jurisdictional issues established therein. The authenticity of this document may be verified at http://www.sgs.com/en/certified-clients-and-products/certified-client-directory. Any unauthorized alteration, forgery or falsification of the content or appearance of this document is unlawful and offenders may be prosecuted to the fullest