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## Note

This image shows a single sheet of white paper with horizontal ruling lines. The lines are evenly spaced and run across the width of the page. There are no margins, text, or other markings on the paper.

**INDEX**

<b>1. Codes .....</b>	<b>pag.18</b>
<b>2. Introduction .....</b>	<b>pag.18</b>
<b>3. Intended use.....</b>	<b>pag.18</b>
<b>4. EC declaration of conformity .....</b>	<b>pag.19</b>
<b>5. General warning .....</b>	<b>pag.19</b>
<b>6. General description.....</b>	<b>pag.19</b>
6.1 Parts view .....	pag.19
<b>7. Assembling/ disassembling .....</b>	<b>pag.20</b>
7.1 Opening.....	pag.20
7.2 Closing .....	pag.21
7.3 Install/remove the armrest.....	pag.21
7.4 Install/remove the swing away and elevating footrest.....	pag.22
7.5 Adjust the swing away legrest .....	pag.22
7.6 Adjust the elevating legrest.....	pag.22
7.7 Install/remove the 30cm wheels .....	pag.23
7.8 Adjust the brakes with 30cm wheels .....	pag.23
7.9 Install/remove the 60cm wheels.....	pag.23
7.10 Adjust the brakes with 60cm wheels.....	pag.24
<b>8. Before use.....</b>	<b>pag.24</b>
<b>9. Using warnings.....</b>	<b>pag.24</b>
<b>10. How to use .....</b>	<b>pag.24</b>
10.1 Get in the wheelchair.....	pag.25
10.2 Get out of the wheelchair.....	pag.25
10.3 Stairs and height difference (how to face them).....	pag.25
10.4 Use the armrest.....	pag.25
10.5 Use the 60cm wheels.....	pag.25
10.6 Use the 30cm wheels.....	pag.25
10.7 Use the casters for narrow ways .....	pag.26
10.8 Carrying .....	pag.26
<b>11. Maintenance .....</b>	<b>pag.26</b>
<b>12. Cleaning and disinfection .....</b>	<b>pag.26</b>
12.1 Cover and stuffing.....	pag.26
12.2 Frame and metal parts.....	pag.26
12.3 Brakes and wheels .....	pag.26
12.4 Disinfection .....	pag.26
<b>13. Conditions of disposal .....</b>	<b>pag.26</b>
<b>14. Spare parts/ accessories .....</b>	<b>pag.26</b>
<b>15. Technical features.....</b>	<b>pag.27</b>
15.1 Dimensions .....	pag.27
15.2 Technical specifications .....	pag.28
<b>16. Warranty.....</b>	<b>pag.29</b>
16.1 Repairing.....	pag.29
16.2 Spare parts.....	pag.29
16.3 Exempt clauses .....	pag.29

**CE** I Class Medical Device  
according to Directive 93/42/EEC  
and further modifications

### 1. Codes

<b>CP100-xx</b>	Series <b>START</b> 60 cm wheels - <b>Self-Propelled</b>
<b>CP102-xx</b>	Series <b>START1</b> 60 cm wheels - <b>Self-Propelled</b>
<b>CP110-xx</b>	Series <b>NEXT</b> 60 cm wheels - <b>Self-Propelled</b>
<b>CP115-xx</b>	Series <b>NEXT GO</b> 30 cm wheels - <b>Transit</b>
<b>CP200-xx</b>	Series <b>DUAL</b> 60 cm wheels - <b>Self-Propelled</b>
<b>CP205-xx</b>	Series <b>DUAL GO</b> 30 cm wheels - <b>Transit</b>
<b>CP300-xx</b>	Series <b>PLUS</b> 60 cm wheels - <b>Self-Propelled - Heavy duty</b>
<b>CP310-xx</b>	Series <b>PLUS GO</b> 30 cm wheels - <b>Transit - Heavy duty</b>
<b>CP500-xx</b>	Series <b>GO!</b> 30 cm wheels - <b>Transit</b>
<b>CP520-xx</b>	Series <b>GO UP!</b> 30 cm wheels - <b>Transit</b>
<b>CP620-xx</b>	Serie <b>SKINNY</b> 50 cm wheels - <b>Small spaces</b>
<b>CP625-xx</b>	Serie <b>SKINNY GO!</b> 30 cm wheels - <b>Transit - Small spaces</b>
<b>CP880-35</b>	Series <b>KIDDY</b> 56 cm wheels - <b>Self-Propelled</b>

(XX =Refers to size available for each wheelchair model ex. 38-40-43-45-46-48-50-55-60 cm)

### 2. Introduction

Thank you for purchasing a **ARDEA ONE** wheelchair by Moretti S.p.A. product range. This user's manual provides some suggestions as how to correctly use the product you have chosen and gives some valuable advice for your safety. Please read through the manual carefully before using the product. Should you have any queries, please contact your retailer for advice and assistance. In case of doubts, we recommend to contact the dealer who can help and advice you properly.

Note: check all parts for shipping damage and test before using. In case of damage, do not use the device. Contact your dealer for further instruction.

### 3. Intended use

This wheelchair is designed for the mobility of people with mobility problems and can be used for self-propelling (if equipped with big rear wheels) or with the aid of an operator.

#### WARNING!

- Do not use the product for a purpose not indicated in this manual
- Moretti S.p.A declines all responsibilities for any consequences resulting from an incorrect use of this product and from unauthorized alteration to the frame of the product
- The manufacturer reserves the right to change the information contained in this document without notice

#### 4. EC declaration of conformity

Moretti S.p.A. hereby declares, under its own exclusive responsibility, that the products manufactured and sold by Moretti S.p.A. and belonging to the wheelchair category are compliant with the applicable provisions outlined in the Medical Devices Directive 93/42/EEC and further modifications. To this end Moretti S.p.A. guarantees and declares as follows, under its own responsibility:

1. The devices in question satisfy the essential requirements of Annex I of **directive 93/42/EEC** as prescribed by Annex VII of the said directive
2. The devices in question **ARE NOT MEASUREMENT DEVICES.**
3. The devices in question **ARE NOT INTENDED FOR CLINICAL INVESTIGATION**
4. The devices in question are sold in **NON STERILE** packing
5. The devices in question should be considered as belonging to class I
6. **DO NOT** use/install the devices in question for any other purpose than that expressly indicated by Moretti S.p.A.
7. Moretti S.p.A. maintains and submits the documentation certifying conformity with **Directive 93/42/EEC** and further modifications to the proper authorities

#### 5. General warning

- **DO NOT** use this product or any available optional equipment without first completely reading and understanding this instruction manual. If you are unable to understand the warnings, cautions or instructions, contact a healthcare professional, dealer or technical personnel before using this equipment otherwise injury or damage may occur.
- Advise with physician or therapist to determinate the correct adjustment and the correct use of the wheelchair
- Keep the packed wheelchair away from heat sources
- **SERVICE LIFE**- the wheelchair use limit is defined by the wear of the parts
- **ALWAYS** keep hands and fingers clear of moving parts to avoid injury
- **DO NOT** allow children to play on or operate the wheelchair
- Close the wheelchair before lifting it.
- **DO NOT** attempt to lift the wheelchair by any removable parts

#### 6. General description

##### 6.1 Parts view

The components can change depending on the model you choose:


A - Frame


B - Wrenches kit for assembling


C - Swing away legrest

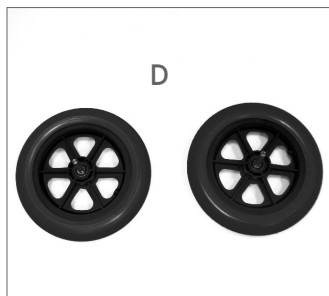
D - 30 cm casters  
(diameter 30 cm) Transit model

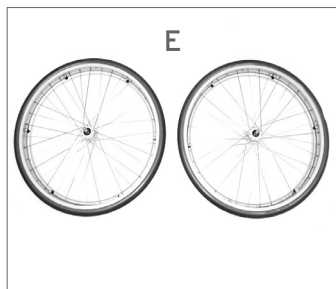E - 60 cm Wheels  
(diameter 60 cm) Self-propelled model

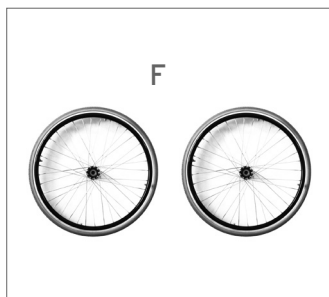
F - 50 cm Wheels for small spaces  
(diameter 50 cm) Self-propelled model

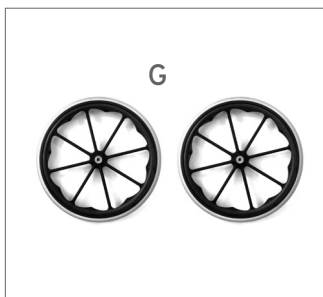G - 56 cm Wheels - Kids model  
(diameter 56 cm) Self-propelled model

## 7. Assembling/disassembling

### 7.1 Opening

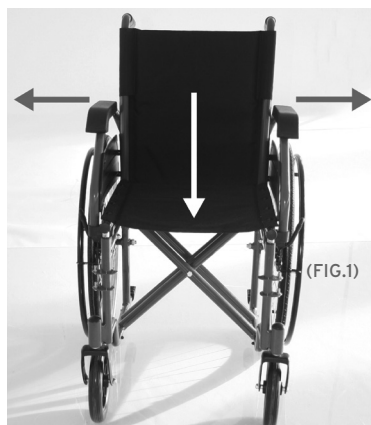Note: make reference to safety cautions in this manual

- Get the handles and pull outward, then press the seat until the wheelchair is completely open (FIG.1)

## 7.2 Closing

- Lift powerfully the seat and close the wheelchair. Push the handles close (FIG.2)

## 7.3 Install/remove the armrest


CP100-xx / CP102-xx / CP110-xx / CP115-xx / CP500-xx / CP520-xx / CP880-35

CP300-xx / CP310-xx / CP200-xx / CP205-xx

CP620-xx / CP625-xx

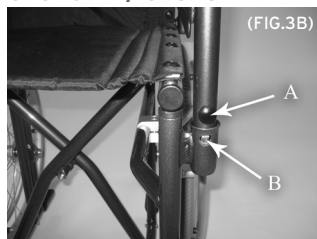A - Push button  
B - Pin for stop


- Place the armrest rear side into the rear, then hook them into the front by overturning them (FIG.3, FIG.3A)
- Insert the tubes in their supports triggering properly locking clip (FIG. 3B)
- To remove, unlock the armrest by pushing the front release lever, then overturn them and extract by pulling the rear release lever (FIG.4, FIG.4A)
- For removal press the release button (A) and then lift to up (FIG. 3B)

#### 7.4 Install/remove the swing away or elevating legrest

- Assembly the swing away or elevating legrest by installing the hinge plates on the footrest onto the hinge pins on the wheelchair frame. Push the legrest toward inside of the wheelchair until it locks into place (FIG.5)
- Remove the legrest by pushing the legrest release lever in ward, rotate outward, lift the legrest off the hinge pins

#### 7.5 Adjust the swing away legrest

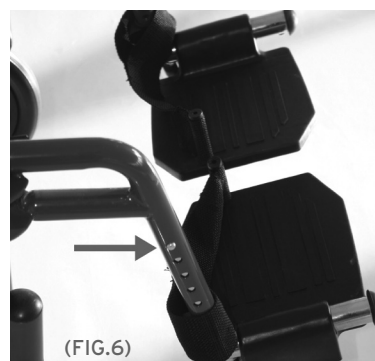- To adjust the length of the footrest, unlock the safety lever at the bottom on the left side, then press the push button, slide the tube until the desired hole, then relock the safety lever. (FIG.6)


#### 7.6 Adjust the elevating legrest (OPTIONAL)

- To adjust the length of the footrest, unlock the safety lever at the bottom on the left side, then press the push button, slide the tube until the desired hole, then relock the safety lever. (FIG.7) (Rif.P)
- Adjust the legrest height by pushing the lever (FIG.7) (Rif.L), and release to lock at the desired height

**7.7 Install/remove the 30cm wheels (only for NEXT, DUAL and SKINNY)**

Place the caster into the hole on the wheelchair frame

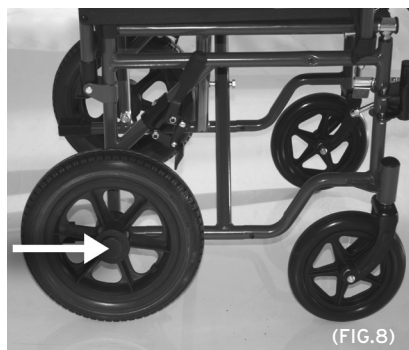
(FIG.8)

**Only for quick release wheels**

- Remove the plastic cap and push the wheel until it locks into the place (FIG.8)
- Remove the wheels: push the button on the centre of the wheel, and extract the wheel (FIG.8)

**Only for wheels with fixed axle**

- Place the small wheel and block it with the nut
- To remove the wheel, unscrew the nut and remove the small wheel

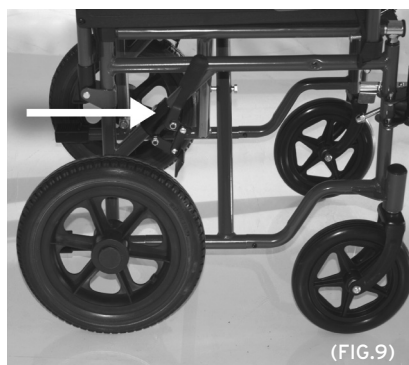**7.8 Adjust the brakes with 30cm wheels (only for NEXT, DUAL and SKINNY)**

(FIG.9)

- When the 30 cm wheels are installed, place the brake on the milling. Lock the device by tightening the Allen wrenches (FIG.9)

**Note adjust the brake device so that it can make the right pressure to securely lock the wheel**

**7.9 Install/remove the 60cm wheels (only for NEXT GO, DUAL GO and SKINNY GO)**

Place the wheel into the central hole on the wheelchair frame

(FIG.10)

**Only for quick release wheels**

- Remove the plastic cap and push the wheel until it locks into the place (FIG.10)
- Remove the wheels: push the button on the centre of the wheel, and extract the wheel (FIG.10)

**Only for wheels with fixed axle (SKINNY GO)**

- Place the big wheel and block it with the nut
- Unscrew the nut and remove the big wheel

### 7.10 Adjust the brakes with 60cm wheels (just for DUAL GO)

- When the 60 cm wheels are installed, place the brake on the milling. Lock the device by tightening the Allen wrenches (FIG.11)

Note adjust the brake device so that it can make the right pressure to securely lock the caster

### 8. Before use

- Inspect the frame for any damage. Check the tubing for cracks and rust, in order to guarantee the safe for people ant thing
- Check the correct assembly, in order to guarantee a safe use
- Check the seat is completely open. If not push the seat down from the lateral frame.
- Check the brake correctly works

### 9. Using warnings

- Check the automotive parts usage, in order to guarantee a safe use of the device for people and things
- Assure that the seat is completely open and locked
- Assure that all parts are correctly assembled. Checking carefully the correct assembling of the wheels
- Assure the brakes correctly works

### 10. How to use

Note: make reference to "GENERAL WARNINGS" chapter

#### **WARNING!**

- Getting in and out of the wheelchair may be performed independently, provided you have a adequate mobility and upper body strength
- Do not stand on removable parts

**DO NOT STAND ON THE FOOTRESTS WHEN GETTING IN OR OUT OF THE WHEELCHAIR, OVERTURNING RISK!!!**

### 10.1 Get in the wheelchair

1. Position the wheelchair as close as possible to you
2. Engage wheel locks . We suggest to remove the armrest to ease the patient getting in the wheelchair
3. Lean on the frame with arms and slowly go down towards the seat

### 10.2 Get out of the wheelchair

1. Position the wheelchair as close as possible along side the seat to which you are transferring
2. Engage wheel locks. We suggest to remove the armrest to ease the patient getting out of the wheelchair
3. Shift body weight into seat with transfer

### 10.3 Stairs and height difference

**WARNING!**

The following methods must be used only from experienced patients. If you are not experienced, let you help from an assistant!

#### GET ON

Position the wheelchair at the curb to be crossed. Balance the wheelchair on the rear wheels, in order to lift the front wheels, enough to lean them on the obstacle. Bend your trunk forward and vigorously push the handrail. Continue pushing until the rear wheels are climbing over the obstacle.

#### GET OFF

To get off the stair follow the reverse procedure, positioning the rear side of the wheelchair at the curb, keep the handrails vigorously and get off the obstacle slowly; then balance the wheelchair on the rear wheels in order to complete the descent

**WARNING!**

We recommend to use this method only in presence of obstacle lower than 10cm!

### 10.4 Use the armrest

The armrest can be flip back, detachable or flip back and detachable

Flip back - this feature is useful when you have to get in or out of the wheelchair.

Push the release lever then turnover the armrest

### 10.5 Use the 60cm wheels

The rear 60cm wheels are used when self-propelled.

Grasp the handrail and wheel through the desired direction

### 10.6 Use the 30cm wheels

The rear 30cm wheels are used when wheel chair is driven by an assistant through the handles

**10.7 Use the narrow casters (only on models with quick release wheels)**

The safety wheels are used to cross narrow spaces (elevator, etc.)

To use the safety wheels remove the rear 60cm wheels

**10.8 Carrying**

The wheelchair has been designed to be easily carried.

Once the armrest and the footrest are removed, can be closed and easily positioned in the car boot, for example. If this can't be possible, position the wheelchair in front of the car rear seat so that it is stable and cannot slip or slide. Lock the wheelchair with the seatbelt

**11. Maintenance**

The **ARDEA ONE** devices by Moretti are checked carefully and supplied with EC mark, once launched on the market. For patient and physician's safety we recommend to let check the device by the manufacturer or by laboratory authorized every 2 years. In case of repair use only original spare parts. Before using follow the "BEFORE USE" chapter.

Moretti **ARDEA ONE** wheelchairs require few elementary and periodic maintenance, below we list the checks to be carried out periodically

braking capacity	monthly
seat cleaning	as needed

**12. Cleaning and disinfection****12.1 Cover and padding**

The cover and padding are in Nylon. Clean with warm water and mild soap.

**12.2 Frame and metallic parts**

Clean with warm water and mild soap. Dry with a soft cloth in order to not scratch the painting.

**12.3 Brakes and wheels**

Clean the wheels with warm water and mild soap. Pay attention in drying the wheels, it's dangerous use the wheelchair with wet wheels.

**12.4 Disinfection**

if you need to disinfect the device use a common disinfectant cleaner

Note: do not use solvent and abrasive product

**13. Disposal**

In case of disposal don't use the inserting container for municipal waste. We recommend to dispose the wheelchair in the appropriate disposal areas for recycle.


**14. Spare parts/accessories**

For spare parts only refer to MORETTI master catalogue

## 15. Technical features

## 15.1 Dimensions

## DIMENSION AND WEIGHT

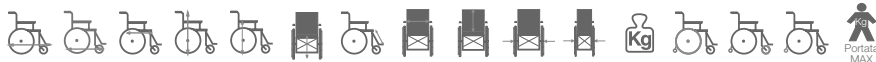CP100-38	113	83	43	90	74	49	58	38	43	56	27	17,5	60	20	-	120
CP100-40								40		58						
CP100-43								43		61						
CP100-45								45		63						
CP100-46								46		63,5						
CP100-48								48		65						
CP102-38	113	83	43	90	74	49	58	38	43	56	27	17,5	60	20	-	120
CP102-40								40		58						
CP102-43								43		61						
CP102-45								45		63						
CP102-46								46		63,5						
CP102-48								48		65						
CP110-40	105	81	43	90	73	49	56	40	42	64	26	18,5	60	20	-	150
CP110-43								43		67						
CP110-45								45		69						
CP110-46								46		69,5						
CP110-48								48		71						
CP110-50								50		74						
CP200-38	113	81	42	90	71	49	59	38	41	58	29	17,5	60	20	4	150
CP200-40								40		60						
CP200-43								43		63						
CP200-45								45		65						
CP200-46								46		65,5						
CP200-48								48		68						
CP200-50	50	70														
CP300-50	118	90	48	92	76	52	59	50	41	73	32	26	60	20	-	200
CP300-55								55		78						
CP300-60								60		83						
CP620-40	100	73	41	92	80	53	47	40	37	51	26	16,5	50	18	-	110
CP620-43								43		54						
CP880-35	94	74	36	86	62	42	38	35	37	56	31	14,5	56	15	-	100

## DIMENSION AND WEIGHT

CP115-40							40		53							
CP115-43							43		56							
CP115-45							45		58							
CP115-46	98	75	41	92	76	50	56	43	58,5	25	18,5	30	20	-	150	
CP115-48									61							
CP115-50									63							
CP205-38								38	52							
CP205-40								40	54							
CP205-43								43	57							
CP205-45	104	76	42	90	71	49	59	45	59	25	17,5	30	20	4	150	
CP205-46								46	59,5							
CP205-48								48	62							
CP205-50								50	64							
CP310-50	104	77	48	92	70	52	59	50	63	30	26	32	20	-	200	
CP310-55								55	68							
CP500-40								40	55							
CP500-43								43	58							
CP500-45	94	70	42	92	72	48	54	45	60	26	17,5	30	20	-	120	
CP500-46								46	60,5							
CP500-48								48	63							
CP520-40								40	55							
CP520-43								43	58							
CP520-45	99	73	42	93	72	48	52	45	60	26	17,5	30	20	-	150	
CP520-46								46	60,5							
CP520-48								48	63							
CP625-40	88	61	41	92	80	53	47	40	51	26	16,5	30	18	-	110	
CP625-43								43	54							

\*The dimensions may change  $\pm 0,5$  cm

## 15.2 Technical specifications

- **Steel frame**
- **Max weight for CP100, CP102, CP500: 120 kg - for CP110, CP115, CP200, CP205, CP520: 150 kg - for CP300, CP310: 200 kg - for CP620, CP625: 110 kg - for CP880-35: 100 kg**
- **Front wheels: 20cm PU**
- **Front wheels (kids model): 15cm PU**
- **Rear small wheels: 30cm PU**
- **Rear big wheels: 60cm PU**
- **Rear big wheels (kids model): 56cm**
- **Fabric: Nylon**
- **Resistance to ignition:** the materials used for cover and padding have been treated with methods specified in Directive EN1021-1/2
- **Reference standards: EN12183**
- **Environmental temperature: -20°C/+ 40°C**
- **Storage:** store in a cool place, away from heat sources, such as sunlight (suggested temperature between +10°C and +40°C)

## 16. Warranty

Moretti products are guaranteed from material or manufacturing faults for 2 years from the purchasing date, except possible exclusion or restriction as follows. The warranty shall not be applied in the possible damages caused by improper use, abuse or alteration, and for the warranty will not be valid if the instructions for use are not strictly adhered to. The correct intended use is specified in this manual. Moretti is not responsible for consequent damages, personal injuries or whatever caused by or in relation to mistaken installation or improper use. Moretti warranty does not cover damages resulting from: natural disaster, not authorized maintenance or repairs, faults caused by problems on electricity supply (when necessary), use of spare parts not covered by Moretti, improper use, alteration not authorized, shipment damages (different from original Moretti shipment), or in case of insufficient maintenance as indicated in the manual. The warranty doesn't cover components subject to wear and tear during the correct use of the device.

### 16.1 Repairing

**Warranty repair:** If a Moretti item presents material or manufacturing faults during the warranty period, Moretti will confirm with customer if the fault can be cover from warranty. Moretti, at its unquestionable discretion, can repair or replace the item, by a Moretti dealer or to Moretti headquarters. Labor cost can be charged to Moretti if the repair is covered by warranty. A repair or a replacement doesn't extend the warranty.

**Repair not cover by warranty:** A product out of warranty can be sent after Moretti authorization. The labor and shipping costs for good out of warranty are to be paid by the customer or by the dealer. The repairs are guaranteed for 6 months from the good received.

**Non defective device:** The customer will be informed if, after the device return and examination, Moretti declares that the device is not faulty. In this case the good will be sent back to customer, the shipping charge will be to be paid by the customer.

### 16.2 Spare parts

The original Moretti spare parts are guaranteed for 6 months from the receiving date.

### 16.3 Exempt clauses

Moretti does not offer any other declaration, explicit or implicit warranty or conditions, including possible declarations, warranties or conditions of merchantability, fitness for a specific purpose, non infringement and non interference, all but what expressly specified in this warranty. Moretti does not guarantee the nonstop and faultless usage. The duration of possible implicit warranties which can be imposed by the law is limited by the warranty period, in the limits of law. Some states or countries don't allow limiting the implicit warranty or the exclusion or the limitation for accidental damages. In such countries, some of those exclusions or limitation may not be applied to the user. The present warranty may be modified without prior notification.

**MORETTI®**

**WARRANTY CERTIFICATE**

**ARDEA®**

**Product** \_\_\_\_\_

**Purchasing date** \_\_\_\_\_

**Authorized dealer** \_\_\_\_\_

**Street** \_\_\_\_\_ **Place** \_\_\_\_\_

**Sold to** \_\_\_\_\_

**Street** \_\_\_\_\_ **Place** \_\_\_\_\_

**MORETTI S.P.A.**

Via Bruxelles, 3 - Melegnano 52022 Cavriglia (Arezzo) Tel. +39 055 96 21 11

[www.morettispa.com](http://www.morettispa.com) email: [info@morettispa.com](mailto:info@morettispa.com)

**MADE IN P.R.C.**

## Notes

This image shows a single page of white paper with horizontal ruling lines. The lines are evenly spaced and run across the width of the page. There are no margins, text, or other markings on the paper.

Via Bruxelles, 3 - Melegnano  
52022 Cavriglia (Arezzo)  
Tel. +39 055 96 21 11

Fax. +39 055 96 21 200  
[www.morettispa.com](http://www.morettispa.com)  
[info@morettispa.com](mailto:info@morettispa.com)