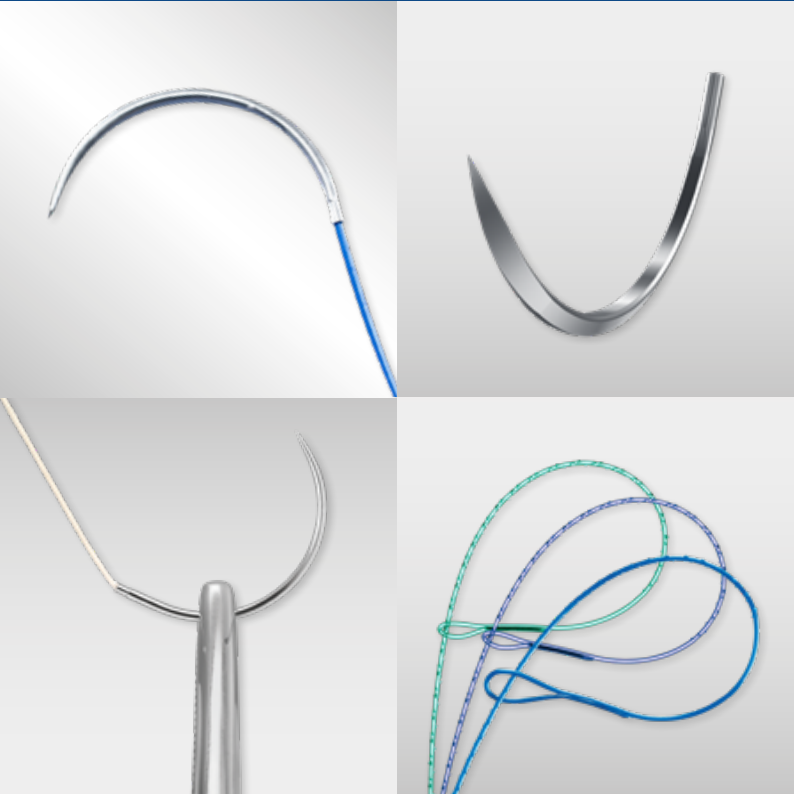




# OPERATING ROOM **ESSENTIALS**

THE COMPLETE  
SOLUTION  
**YOU CAN TRUST**

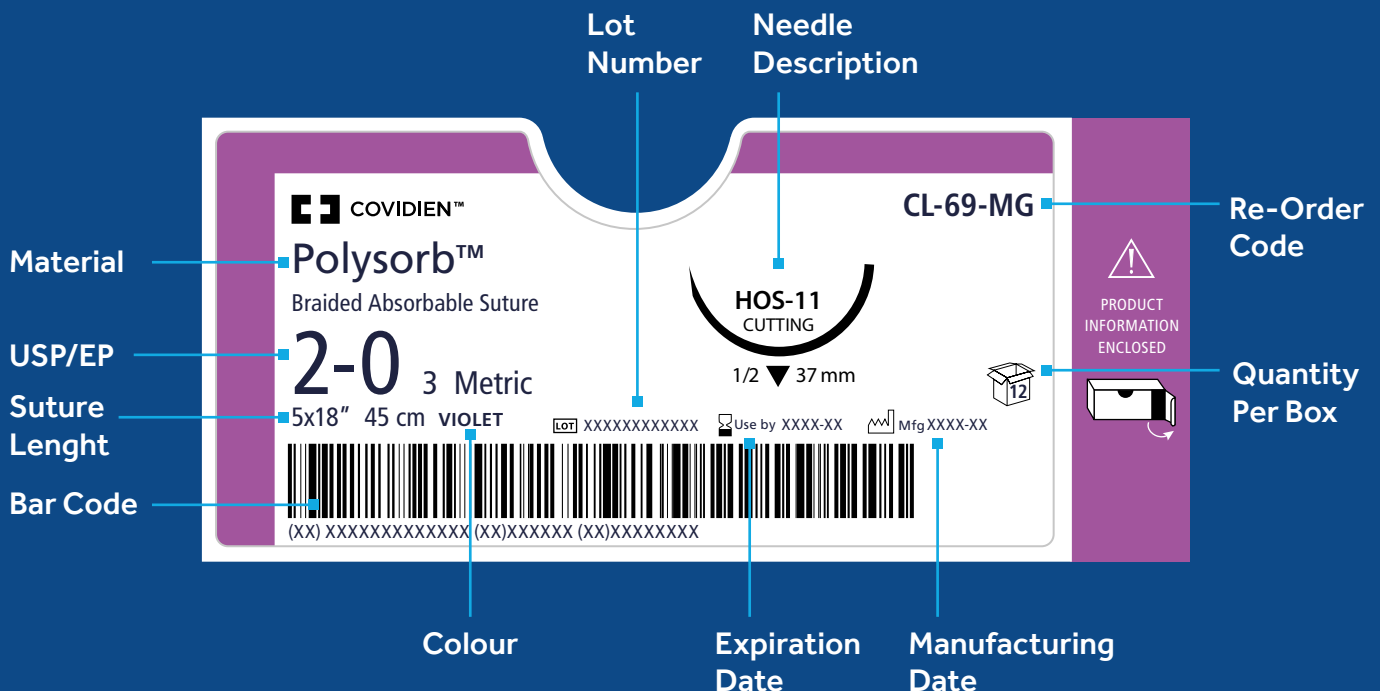
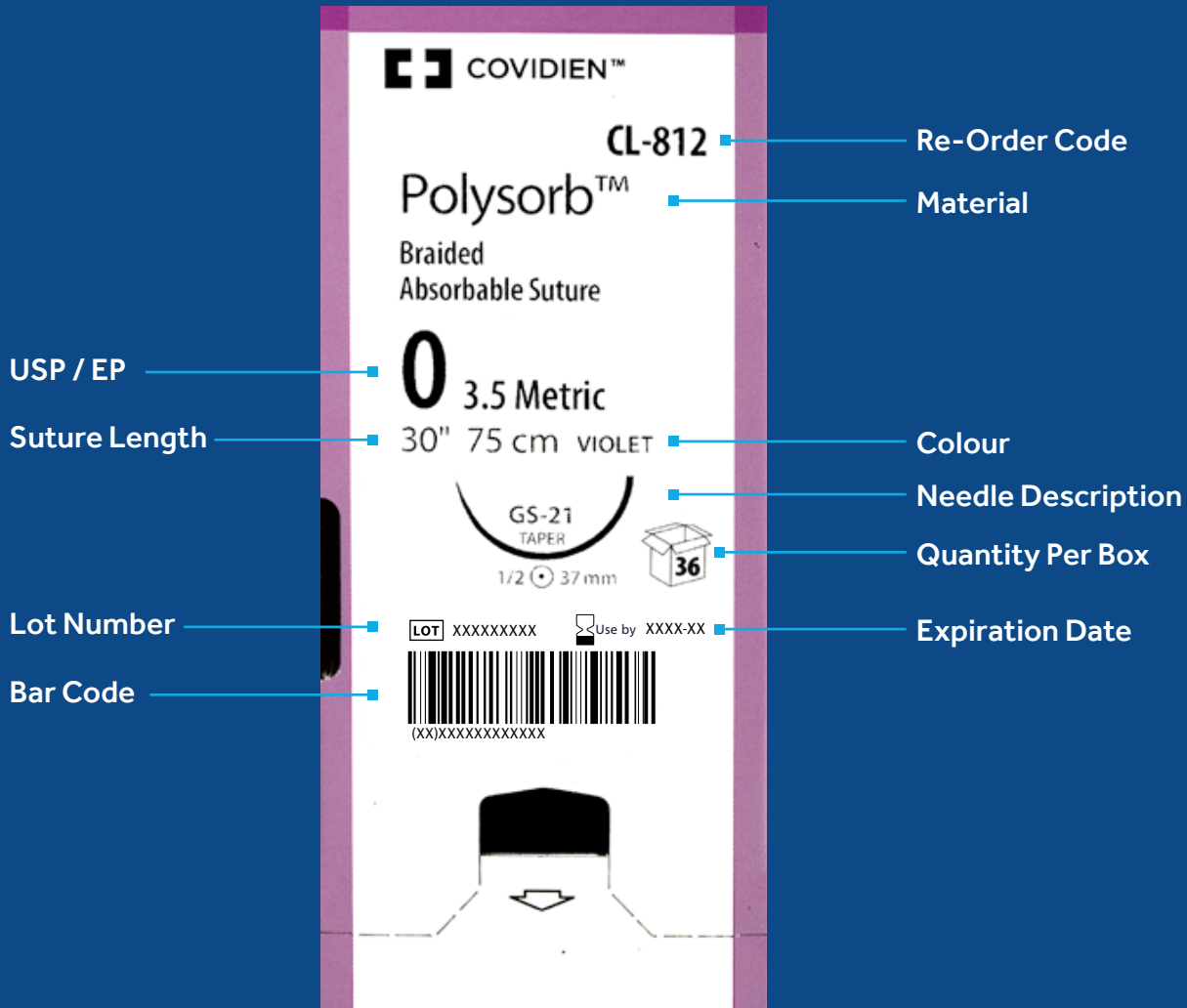


Wound Closure  
**Catalogue**

October 2019

**Medtronic**  
Further, Together

# HOW TO READ A SUTURE LABEL



# WHY MEDTRONIC?

Medtronic have an extensive heritage in partnerships with Public Hospitals and Private Groups. Delivering value-based healthcare, clinical expertise and innovation through products and services that are focused on patient care and clinical outcomes.

## OR ESSENTIALS

We pride ourselves on being a brand you can trust. With our proven track record of **quality products and services**, you can feel reassured that Medtronic's OR Essentials will help you to continue providing **the very best patient care** from the first mark to the last stitch. Let's embrace the future together.



Laparoscopic Access



Surgical Instruments



Wound Closure

# INTRODUCTION

## A Comprehensive Portfolio of Wound Closure Materials and Accessories

We're all part of the same community. A community entrusted with enhancing and saving the lives of patients. At Medtronic, we actively partner with you to provide your patients with the most optimal care. Every day. Whatever role you play, we're here for you — leveraging over 100 years of expertise and insight to innovate wound closure solutions. We are committed to remaining by your side for the next 100 years, evolving with your changing needs so you can keep delivering the highest levels of patient care. Our comprehensive portfolio of effective wound closure solutions is for the full spectrum of general and specialty procedures.

### Absorbable Sutures

Absorbable sutures from Medtronic provide strength retention and absorption profiles to meet your surgical specialty and tissue healing requirements.

### Non-Absorbable Sutures

Non-absorbable sutures from Medtronic can be used in a range of procedures, including cardiovascular, ophthalmic and neurological surgery.





















### Barbed Wound Closure

The V-Loc™ wound closure device is a revolutionary technology that eliminates the need to tie knots, so you can close incisions up to 50% faster without compromising strength and security

### Wound Closure Accessories

Wound closure accessories from Medtronic include bone wax, bolsters, suture racks, cardiovascular pledgets and umbilical tape

# CONTENTS

	Velosorb™ Fast .....	6
	Polysorb™ .....	18
	Caprosyn™ .....	50
	Biosyn™ .....	60
	Maxon™ / Maxon™ CV .....	76
	Surgipro™ / Surgipro™ II .....	94
	Ti-Cron™ / Surgidac™ .....	114
	Vascufil™ .....	136
	Flexon™ .....	140
	Steel .....	144
	Dermalon™ / Monosof™ .....	150
	Novafil™ .....	168
	Surgilon™ .....	178
	Sofsilk™ .....	186
	V-Loc™ 90 Absorbable .....	204
	V-Loc™ 180 Absorbable .....	210
	V-Loc™ PBT .....	216
	Pledgets .....	220
	Wound closure accessories .....	223
	Needles chart .....	224
	Ordering information .....	228

# VELOSO RB™ FAST

## Braided Fast Absorbable Suture

There when you need it. Gone when you don't.





Please refer to the individual package IFUs for complete instructions, indications, contraindications, warnings and precautions.

<b>Description</b>	Velosorb™ Fast absorbable braided sutures are composed of an absorbable synthetic polyester comprised of glycolide and lactide (derived from glycolic and lactic acids). Velosorb™ Fast absorbable sutures are prepared by coating the suture with a mixture of a glycolide and lactide copolymer and calcium stearate.
<b>Indications</b>	Velosorb™ Fast absorbable sutures are indicated for use in soft tissue approximation of the skin, mucosa and conjunctiva, where only 7-10 days of wound support is required, but not for use in ligation, cardiovascular or neurological tissue.
<b>Actions</b>	Velosorb™ Fast absorbable sutures elicit a minimal acute inflammatory reaction in tissue, which is followed by a gradual encapsulation of the suture by fibrous connective tissue. Progressive loss of tensile strength and eventual absorption of Velosorb™ Fast absorbable sutures occurs by means of hydrolysis, where the glycolide/lactide copolymer is broken down to glycolic and lactic acids which are subsequently absorbed and metabolized by the body.
<b>Tensile Strength</b>	Studies indicate tensile strength averages for Velosorb™ Fast absorbable sutures are no less than 45% of its original tensile strength 5 days post implantation. All of the original tensile strength is essentially gone by 14 days.
<b>Sterilization</b>	Gamma irradiation.
<b>Absorbtion</b>	Absorption begins as a loss of tensile strength without appreciable loss of mass. Absorption of Velosorb™ Fast absorbable sutures is essentially complete between the 40 <sup>th</sup> and 50 <sup>th</sup> day.
<b>Packaging</b>	Velosorb™ Fast absorbable sutures are available in USP sizes 1 (4 Metric) through 6-0 (0.7 Metric). They are available undyed (natural). The sutures are supplied sterile, in pre-cut lengths affixed to various needle types. Please refer to the individual package IFUs for complete instructions, indications, contraindications, warnings and precautions.

## Velosorb™ Fast absorbable suture






Fast braided absorbable suture for use in soft tissue approximation of the skin and mucosa, when only 7–10 days of wound support is required.

PRE-CUT							
Needles Specifications	Size		Length cm	Color	Units Per box	Comments	Reorder code
	USP	EP					
	<b>4-0</b>	1.5	6 × 45	Undyed	36		CV-8977
	<b>4-0</b>	1.5	3 × 45	Undyed	36		CV-8643
	<b>3-0</b>	2	6 × 45	Undyed	36		CV-8978
	<b>3-0</b>	2	3 × 45	Undyed	36		CV-8644
	<b>2-0</b>	3	6 × 45	Undyed	36		CV-8979
	<b>2-0</b>	3	5 × 75	Undyed	36		CV-8653
	<b>0</b>	3.5	6 × 45	Undyed	36		CV-8980






REVERSE CUTTING							
Needles Specifications	Size		Length cm	Color	Units Per box	Comments	Reorder code
	USP	EP					
<b>C-1</b> 3/8 Circle 12 mm 	<b>5-0</b>	1	75	Undyed	36		SV-2289
	<b>4-0</b>	1.5	45	Undyed	36		SV-1290
	<b>4-0</b>	1.5	75	Undyed	36		SV-2290
<b>C-12</b> 3/8 Circle 16 mm 	<b>5-0</b>	1	75	Undyed	36		SV-2296
	<b>4-0</b>	1.5	75	Undyed	36		SV-2294
	<b>3-0</b>	2	75	Undyed	36		SV-2295










## REVERSE CUTTING

Needles Specifications	Size		Length cm	Color	Units Per box	Comments	Reorder code
	USP	EP					
<b>C-13</b> 3/8 Circle 19 mm 	<b>5-0</b>	1	45	Undyed	36		SV-2910
	<b>4-0</b>	1.5	45	Undyed	36		SV-2920
	<b>4-0</b>	1.5	75	Undyed	36		SV-2297
	<b>3-0</b>	2	45	Undyed	36		SV-2930
	<b>3-0</b>	2	75	Undyed	36		SV-2298
	<b>2-0</b>	3	75	Undyed	36		SV-2299
<b>C-14</b> 3/8 Circle 24 mm 	<b>4-0</b>	1.5	75	Undyed	36		SV-2251
	<b>3-0</b>	2	45	Undyed	36		SV-4431
	<b>3-0</b>	2	75	Undyed	36		SV-2252
	<b>2-0</b>	3	75	Undyed	36		SV-2253
<b>C-15</b> 3/8 Circle 26 mm 	<b>3-0</b>	2	75	Undyed	12		SV-9940G
	<b>2-0</b>	3	75	Undyed	12		SV-9941G
<b>C-16</b> 3/8 Circle 30 mm 	<b>3-0</b>	2	90	Undyed	36		SV-2285
	<b>2-0</b>	3	90	Undyed	36		SV-2286
	<b>0</b>	3.5	90	Undyed	36		SV-2287
<b>C-17</b> 3/8 Circle 39 mm 	<b>2-0</b>	3	75	Undyed	12		SV-9961G


## REVERSE CUTTING

Needles Specifications	Size		Length cm	Color	Units Per box	Comments	Reorder code
	USP	EP					
<b>C-22</b> 1/2 Circle 18 mm 	4-0	1.5	45	Undyed	36		SV-4550
<b>C-23</b> 1/2 Circle 24 mm 	4-0	1.5	45	Undyed	36		SV-9926
	3-0	2	45	Undyed	36		SV-9927
	2-0	3	45	Undyed	36		SV-9928
<b>GS-10</b> 1/2 Circle 26 mm 	2-0	3	75	Undyed	12		CV-9936G
<b>GS-11</b> 1/2 Circle 37 mm 	3-0	2	135	Undyed	12		CV-931G
	2-0	3	90	Undyed	36		CV-4721
	2-0	3	135	Undyed	12		CV-932G
	0	3.5	90	Undyed	36		CV-4731
<b>GS-12</b> 1/2 Circle 40 mm 	0	3.5	90	Undyed	36		CV-485
	1	4	90	Undyed	36		CV-486




## PREMIUM REVERSE CUTTING

Needles Specifications	Size		Length cm	Color	Units Per box	Comments	Reorder code
	USP	EP					
<b>P-10</b> 3/8 Circle 11 mm 	6-0	0.7	45	Undyed	36		SV-9913
	6-0	0.7	45	Violet	36		SV-32
	6-0	0.7	75	Undyed	36		SV-4893
	5-0	1	45	Undyed	36		SV-9915
<b>P-13</b> 3/8 Circle 13 mm 	5-0	1	45	Undyed	36		SV-493
	4-0	1.5	45	Undyed	36		SV-494
	4-0	1.5	75	Undyed	36	Surgalloy™	SVP-494
<b>P-11</b> 3/8 Circle 16 mm 	5-0	1	75	Undyed	36	Surgalloy™	SVP-2296
	4-0	1.5	75	Undyed	36	Surgalloy™	SVP-2294
	3-0	2	75	Undyed	36	Surgalloy™	SVP-2295
<b>P-12</b> 3/8 Circle 19 mm 	4-0	1.5	45	Undyed	12		SV-496G
	4-0	1.5	45	Undyed	36		SV-496
	4-0	1.5	75	Undyed	36		SV-426
	4-0	1.5	90	Undyed	36		SV-412
	3-0	2	45	Undyed	12		SV-497G
	3-0	2	75	Undyed	36		SV-425
	2-0	3	75	Undyed	36		SV-424
<b>P-14</b> 3/8 Circle 24 mm 	4-0	1.5	75	Undyed	36		SVP-2251
	3-0	2	45	Undyed	36		SV-6830
	3-0	2	75	Undyed	36		SV-935
	2-0	3	45	Undyed	36		SV-6840
	2-0	3	75	Undyed	36		SV-439
<b>P-17</b> 3/8 Circle 30 mm 	3-0	2	90	Undyed	36		SVP-2285
<b>P-18</b> 3/8 Circle 38 mm 	3-0	2	75	Undyed	12		SV-9937G
	2-0	3	75	Undyed	12		SV-9938G





## PREMIUM REVERSE CUTTING

Needles Specifications	Size		Length cm	Color	Units Per box	Comments	Reorder code
	USP	EP					
<b>P-24</b> 1/2 Circle    16 mm 	<b>5-0</b>	1	45	Undyed	12		SV-9910G
	<b>4-0</b>	1.5	45	Undyed	12		SV-9911G







## CONVENTIONAL CUTTING

Needles Specifications	Size		Length cm	Color	Units Per box	Comments	Reorder code
	USP	EP					
<b>PC-10</b> 3/8 Circle    13 mm 	<b>5-0</b>	1	45	Undyed	12		SV-834G
	<b>4-0</b>	1.5	45	Undyed	12		SV-835G
<b>PC-11</b> 3/8 Circle    16 mm 	<b>6-0</b>	0.7	75	Undyed	12		SV-489G
	<b>5-0</b>	1	45	Undyed	12		SV-844G
	<b>4-0</b>	1.5	45	Undyed	12		SV-845G
	<b>4-0</b>	1.5	75	Undyed	12		SV-488G
	<b>3-0</b>	2	75	Undyed	12		SV-487G
<b>SC-11</b> Straight    11 mm 	<b>4-0</b>	1.5	45	Undyed	12		SV-9909G


## CONVENTIONAL CUTTING

Needles Specifications	Size		Length cm	Color	Units Per box	Comments	Reorder code
	USP	EP					
<b>SC-1</b> Straight    51 mm  	<b>4-0</b>	1.5	75	Undyed	36		SV-2238
	<b>3-0</b>	2	75	Undyed	36		SV-2239
	<b>2-0</b>	3	75	Undyed	36		SV-2240
<b>SC-2</b> Straight    60 mm  	<b>3-0</b>	2	75	Undyed	36		SV-9943
	<b>2-0</b>	3	75	Undyed	12		SV-9944G

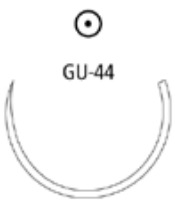
## TAPER POINT

Needles Specifications	Size		Length cm	Color	Units Per box	Comments	Reorder code
	USP	EP					
<b>CV-1</b> 3/8 Circle    9 mm  	<b>8-0</b>	0.4	30	Violet	12		CV-9916G
<b>CV-23</b> 1/2 Circle    17 mm  	<b>5-0</b>	1	75	Undyed	36		CV-1213
	<b>4-0</b>	1.5	75	Undyed	36		CV-1214
	<b>3-0</b>	2	75	Undyed	36		CV-1215
<b>CV-25</b> 1/2 Circle    22 mm  	<b>4-0</b>	1.5	75	Undyed	36	Optivis™	CV-3100
	<b>3-0</b>	2	75	Undyed	36	Optivis™	CV-3110
	<b>2-0</b>	3	75	Undyed	36		CV-3120



## TAPER POINT

Needles Specifications	Size		Length cm	Color	Units Per box	Comments	Reorder code
	USP	EP					
<b>V-20</b> 1/2 Circle    26 mm 	<b>4-0</b>	1.5	75	Undyed	36		CV-415
	<b>3-0</b>	2	75	Undyed	12		CV-416G
	<b>2-0</b>	3	75	Undyed	12		CV-417G
	<b>0</b>	3.5	75	Undyed	36		CV-418
<b>GS-22</b> 1/2 Circle    27 mm 	<b>2-0</b>	3	75	Undyed	36		CV-55
	<b>0</b>	3.5	75	Undyed	36		CV-56
	<b>1</b>	4	75	Undyed	36		CV-57
<b>V-30</b> 1/2 Circle    30 mm 	<b>0</b>	3.5	75	Undyed	36		CV-318
	<b>1</b>	4	75	Undyed	36		CV-319
<b>GS-20</b> 1/2 Circle    34 mm 	<b>2-0</b>	3	75	Undyed	12		CV-9975G
<b>GS-21</b> 1/2 Circle    37 mm 	<b>3-0</b>	2	75	Undyed	36		CV-322
	<b>3-0</b>	2	90	Undyed	36		CV-944
	<b>2-0</b>	3	75	Undyed	36		CV-323
	<b>2-0</b>	3	90	Undyed	36		CV-945
	<b>0</b>	3.5	75	Undyed	36		CV-1324
	<b>0</b>	3.5	90	Undyed	36		CV-946
	<b>1</b>	4	90	Undyed	36		CV-947

## TAPER POINT


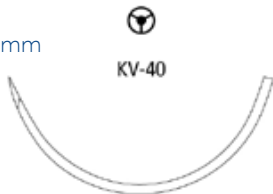
Needles Specifications	Size		Length cm	Color	Units Per box	Comments	Reorder code
	USP	EP					
<b>GS-24</b> 1/2 Circle    40 mm 	<b>2-0</b>	3	90	Undyed	36		CV-9592
	<b>0</b>	3.5	90	Undyed	36		CV-9591
	<b>1</b>	4	90	Undyed	36		CV-9590
<b>GS-25</b> 1/2 Circle    48 mm 	<b>2-0</b>	3	90	Undyed	36		CV-606
	<b>0</b>	3.5	90	Undyed	36		CV-604
	<b>1</b>	4	90	Undyed	36		CV-605
<b>GU-45</b> 5/8 Circle    37 mm 	<b>2-0</b>	3	75	Undyed	36		UV-412
	<b>1</b>	4	90	Undyed	36		UV-401
<b>GU-44</b> 5/8 Circle    40 mm 	<b>0</b>	3.5	90	Undyed	36		UV-1381
	<b>1</b>	4	90	Undyed	36		UV-1382
<b>ST-4</b> Straight    19 mm 	<b>4-0</b>	1.5	75	Undyed	36		CV-1673

## TAPERCUTTING


Needles Specifications	Size USP	EP	Length cm	Color	Units Per box	Comments	Reorder code
<b>KV-11</b> 3/8 Circle 13 mm 	5-0	1	75	Undyed	36		CV-2293
<b>KV-15</b> 3/8 Circle 17 mm 	4-0	1.5	45	Undyed	12		CV-2275G
	4-0	1.5	75	Undyed	36		CV-2279
	3-0	2	75	Undyed	36		CV-2280
<b>KV-5</b> 1/2 Circle 17 mm 	3-0	2	75	Undyed	12		CV-1674G
<b>KV-7</b> 1/2 Circle 26 mm 	3-0	2	75	Undyed	12		CV-2003G
	2-0	3	90	Undyed	36		CV-1002
<b>KV-34</b> 1/2 Circle 37 mm 	3-0	2	90	Undyed	36		CV-916
	2-0	3	90	Undyed	36		CV-917
	2-0	3	120	Undyed	36		CV-927
	0	3.5	90	Undyed	36		CV-919
	1	4	90	Undyed	36		CV-918
<b>KV34-C14</b> 1/2 Circle 40 mm Double-arm	2-0	3.5/3	2×90	Undyed	24	2 strands KV-34	CV-0146
	0	3.5/3	1×70	Undyed	24	1 strand C-14	Episiotomy set




## TAPERCUTTING

Needles Specifications	Size		Length cm	Color	Units Per box	Comments	Reorder code
	USP	EP					
<b>KV-37</b> 1/2 Circle 40 mm 	<b>2-0</b>	3	90	Undyed	36		CV-2218
	<b>0</b>	3.5	90	Undyed	36		CV-2217
	<b>1</b>	4	90	Undyed	36		CV-2216
<b>KV37-C-14</b> 1/2 Circle 40.0 mm Double-arm	<b>2-0</b>	3.5/3	1 x 70	Undyed	24	1 strand C-14	CV-0046
	<b>0</b>	3.5/3	2 x 90	Undyed	24	2 strands KV-37	Episiotomy set
<b>KV-40</b> 1/2 Circle 48 mm 	<b>2-0</b>	3	90	Undyed	36		CV-7192
	<b>0</b>	3.5	90	Undyed	36		CV-7191
	<b>1</b>	4	90	Undyed	36		CV-7190

## BLUNT POINT

Needles Specifications	Size		Length cm	Color	Units Per box	Comments	Reorder code
	USP	EP					
<b>BGS-24</b> 1/2 Circle 40 mm 	<b>0</b>	3.5	90	Undyed	36		CV-9953

## SPATULA

Needles Specifications	Size		Length cm	Color	Units Per box	Comments	Reorder code
	USP	EP					
<b>SE-140-8</b> 3/8 Circle 6.3 mm Curvature 140° 	<b>8-0</b>	0.4	30	Violet	12		CV-501G
	<b>7-0</b>	0.5	30	Violet	12		CV-500G

# POLYSORB™

## Braided Absorbable Suture

Stronger<sup>1</sup> during the critical wound healing period



1. van Ramshorst GH et al. World J Surg 2010;54:20–7

Please refer to the individual package IFUs for complete instructions, indications, contraindications, warnings and precautions.

<b>Description</b>	Polysorb™ braided sutures are composed of Lactomer™ glycolide/lactide copolymer which is a synthetic polyester composed of glycolide and lactide (derived from glycolic and lactic acids). Polysorb™ sutures are prepared by coating the suture with a mixture of a caprolactone/glycolide copolymer and calcium stearoyl lactylate.
<b>Indications</b>	Polysorb™ sutures are indicated for use in soft tissue approximation or ligation and ophthalmic surgery, but not in cardiovascular or neural tissue.
<b>Actions</b>	Polysorb™ sutures elicit a minimal acute inflammatory reaction in tissue, which is followed by a gradual encapsulation of the suture by fibrous connective tissue. Progressive loss of tensile strength and eventual absorption of Polysorb™ sutures occurs by means of hydrolysis, where the Lactomer™ glycolide/lactide copolymer is broken down to glycolic and lactic acids which are subsequently absorbed and metabolized by the body.
<b>Tensile Strength</b>	Studies indicate tensile strength averages for Polysorb™ sutures are approximately 140% of USP and E.P. minimum knot strength initially, are approximately 80% after two weeks and in excess of 30% at three weeks post implant.
<b>Sterilization</b>	Ethylene oxide.
<b>Absorbtion</b>	Absorption begins as a loss of tensile strength without appreciable loss of mass. Absorption of Polysorb™ sutures is essentially complete between the 56 <sup>th</sup> and 70 <sup>th</sup> day.
<b>Packaging</b>	Polysorb™ sutures are available in USP sizes 2 (5 Metric) through 8-0 (0.4 Metric). They are available undyed (natural) or violet colored. The sutures are supplied sterile, in pre-cut lengths and ligating reels, non-neededled or affixed to various needle types using both permanent and removable needle attachment techniques. Please refer to the individual package IFUs for complete instructions, indications, contraindications, warnings and precautions.

# Polysorb™

Polysorb™ suture is a brand you have trusted for over 20 years with a history of improvements and performance.

PRE-CUT							
Needles Specifications	Size		Length cm	Color	Units Per box	Comments	Reorder code
	USP	EP					
	4-0	1.5	3×45	Violet	36	Pre-cut	L-71
	4-0	1.5	6×45	Violet	24	Pre-cut	L-20
	4-0	1.5	12×45	Violet	24	Pre-cut	L-11
	4-0	1.5	12×45	Undyed	24	Pre-cut	L-1-
	3-0	2	3×45	Violet	36	Pre-cut	L-72
	3-0	2	6×45	Violet	24	Pre-cut	L-22
	3-0	2	12×45	Undyed	24	Pre-cut	L-2-
	3-0	2	12×45	Violet	24	Pre-cut	L-12
	3-0	2	12×60	Violet	24	Pre-cut	L-32
	3-0	2	3×75	Violet	24	Pre-cut	L-402
	3-0	2	6×75	Violet	24	Pre-cut	L-38
	2-0	3	3×45	Violet	36	Pre-cut	L-73
	2-0	3	12×45	Undyed	24	Pre-cut	L-3
	2-0	3	6×45	Violet	24	Pre-cut	L-23
	2-0	3	12×45	Violet	24	Pre-cut	L-13-
	2-0	3	12×60	Violet	24	Pre-cut	L-33
	2-0	3	12×75	Undyed	24	Pre-cut	L-43
	2-0	3	3×75	Violet	24	Pre-cut	L-403
	2-0	3	6×75	Violet	24	Pre-cut	L-39
	0	3.5	3×45	Violet	36	Pre-cut	L-74
	0	3.5	6×45	Undyed	24	Pre-cut	L-4
	0	3.5	12×45	Violet	24	Pre-cut	L-27
	0	3.5	12×45	Undyed	24	Pre-cut	L-24
	0	3.5	6×45	Violet	24	Pre-cut	L-14
	0	3.5	12×60	Violet	24	Pre-cut	L-34
	0	3.5	3×75	Violet	24	Pre-cut	L-404
	0	3.5	6×75	Violet	24	Pre-cut	L-40
	1	4	3×45	Violet	36	Pre-cut	L-75
	1	4	6×45	Undyed	24	Pre-cut	L-5
	1	4	12×45	Violet	24	Pre-cut	L-35
	1	4	12×45	Undyed	24	Pre-cut	L-25
	1	4	6×45	Violet	24	Pre-cut	L-15
	1	4	12×60	Violet	24	Pre-cut	L-36

## PRE-CUT

Needles Specifications	Size		Length cm	Color	Units Per box	Comments	Reorder code
	USP	EP					
	<b>1</b>	4	6×75	Violet	24	Pre-cut	L-41
	<b>1</b>	4	3×75	Violet	24	Pre-cut	L-405
	<b>2</b>	5	3×45	Violet	36	Pre-cut	L-76
	<b>2</b>	5	12×45	Violet	24	Pre-cut	L-28
	<b>2</b>	5	12×45	Undyed	24	Pre-cut	L-26
	<b>2</b>	5	6×75	Violet	24	Pre-cut	L-42
	<b>2</b>	5	3×75	Violet	24	Pre-cut	L-406

## STANDARD LENGTH

Needles Specifications	Size		Length cm	Color	Units Per box	Comments	Reorder code
	USP	EP					
	<b>4-0</b>	1.5	150	Violet	36	Standard length	L-111
	<b>3-0</b>	2	150	Undyed	36	Standard length	L-102
	<b>3-0</b>	2	150	Violet	36	Standard length	L-112
	<b>2-0</b>	3	150	Undyed	36	Standard length	L-103
	<b>2-0</b>	3	150	Violet	36	Standard length	L-113
	<b>0</b>	3.5	150	Violet	36	Standard length	L-114
	<b>0</b>	3.5	150	Undyed	36	Standard length	L-104
	<b>2</b>	5	150	Violet	36	Standard length	L-116
	<b>1</b>	4	150	Violet	36	Standard length	L-115





## REEL

Needles Specifications	Size		Length cm	Color	Units Per box	Comments	Reorder code
	USP	EP					
	4-0	1.5	150	Violet	24	Reel	LL-111
	4-0	1.5	150	Undyed	24	Reel	LL-101
	4-0	1.5	250	Violet	24	Reel	LL-221
	3-0	2	150	Violet	24	Reel	LL-112
	3-0	2	150	Undyed	24	Reel	LL-102
	3-0	2	250	Undyed	24	Reel	LL-322
	3-0	2	250	Violet	24	Reel	LL-222
	2-0	3	150	Violet	24	Reel	LL-113
	2-0	3	150	Undyed	24	Reel	LL-103
	2-0	3	250	Undyed	24	Reel	LL-323
	2-0	3	250	Violet	24	Reel	LL-223
	2-0	3	366	Violet	24	Reel	LL-233
	0	3.5	150	Violet	24	Reel	LL-114
	0	3.5	150	Undyed	24	Reel	LL-104
	0	3.5	250	Violet	24	Reel	LL-224
	1	4	150	Violet	24	Reel	LL-115
	1	4	250	Violet	24	Reel	LL-225
	2	5	150	Violet	24	Reel	LL-116
	2	5	250	Violet	24	Reel	LL-226

## REVERSE CUTTING


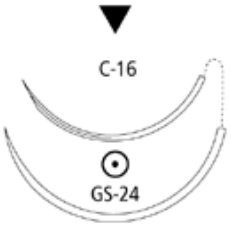



Needles Specifications	Size USP	EP	Length cm	Color	Units Per box	Comments	Reorder code
<b>HE-6</b> 3/8 Circle 7.4 mm Double-arm Curvature 110°	6-0	0.7	45	Violet	12	Surgalloy™	L-1551-K
<b>C-1</b> 3/8 Circle 12 mm	6-0	0.7	45	Undyed	36		SL-171
		5-0	1	45	Undyed	36	
<b>HE-3</b> 3/8 Circle 13 mm Double-arm Curvature 155°	5-0	1	30	Undyed	12	Surgalloy™	L-4323K

## REVERSE CUTTING




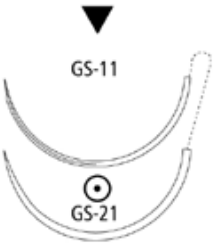
Needles Specifications	Size USP	EP	Length cm	Color	Units Per box	Comments	Reorder code
<b>C-12</b> 3/8 Circle 16 mm 	5-0	1	75	Violet	36		SL-630
	4-0	1.5	75	Violet	36		SL-607
	4-0	1.5	75	Undyed	36		SL-605
	3-0	2	75	Violet	36		SL-608
<b>C-13</b> 3/8 Circle 19 mm 	5-0	1	45	Undyed	36		SL-639
	5-0	1	75	Violet	36		SL-634
	5-0	1	75	Undyed	36		SL-690
	4-0	1.5	45	Undyed	12		SL-638
	4-0	1.5	75	Undyed	12		SL-691G
	4-0	1.5	75	Undyed	36		SL-691
	4-0	1.5	75	Violet	12		SL-635G
	4-0	1.5	75	Violet	36		SL-635
	3-0	2	45	Undyed	36		SL-637
	3-0	2	75	Undyed	36		SL-822
	3-0	2	75	Undyed	12		SL-822G
	3-0	2	75	Violet	36		SL-636
	<b>C-14</b> 3/8 Circle 24 mm 	4-0	1.5	75	Violet	36	
4-0		1.5	75	Undyed	36		SL-692
3-0		2	75	Undyed	36		SL-693
3-0		2	75	Violet	36		SL-694
2-0		3	75	Undyed	36		SL-823
2-0		3	75	Violet	36		CL-823
<b>C-15</b> 3/8 Circle 26 mm 	3-0	2	5 × 45	Undyed	12	D-Tach™	SL-73-MG
	2-0	3	75	Undyed	36		SL-586
	0	3.5	75	Undyed	36		SL-587






## REVERSE CUTTING

Needles Specifications	Size		Length cm	Color	Units Per box	Comments	Reorder code
	USP	EP					
<b>C-16</b> 3/8 Circle 30 mm 	<b>3-0</b>	2	75	Undyed	36		SL-696
	<b>2-0</b>	3	75	Violet	36		SL-714
	<b>2-0</b>	3	75	Undyed	36		SL-611
	<b>0</b>	3.5	75	Violet	36		SL-713
	<b>0</b>	3.5	75	Undyed	36		SL-610
	<b>0</b>	3.5	90	Violet	36		SL-712
<b>C-16 - GS-24</b> 3/8 Circle 30 mm Double-arm (Not to scale) 	<b>3-0</b>	2	150	Undyed	36		CL-976
	<b>2-0</b>	3	150	Undyed	36		CL-977
<b>C-17</b> 3/8 Circle 39 mm 	<b>2-0</b>	3	75	Undyed	36		SL-674
	<b>2</b>	5	90	Undyed	36		SL-675
<b>C-23</b> 1/2 Circle 24 mm 	<b>3-0</b>	2	45	Undyed	36		CL-85
	<b>3-0</b>	2	5 × 45	Undyed	12	D-Tach™	SL-85-MG
	<b>3-0</b>	2	5 × 45	Violet	12	D-Tach™	SL-83-MG
	<b>3-0</b>	2	75	Violet	36		SL-632
	<b>3-0</b>	2	75	Undyed	36		SL-613
	<b>2-0</b>	3	75	Violet	36		SL-631
	<b>2-0</b>	3	75	Undyed	36		SL-612
<b>GS-10</b> 1/2 Circle 26 mm 	<b>3-0</b>	2	5 × 45	Undyed	12	D-Tach™	CL-53MG
	<b>3-0</b>	2	75	Undyed	36		CL-868
	<b>2-0</b>	3	5 × 45	Undyed	12	D-Tach™	CL-52-MG
	<b>2-0</b>	3	75	Undyed	36		CL-869
	<b>2-0</b>	3	75	Violet	36		CL-769
	<b>0</b>	3.5	75	Undyed	36		CL-870
	<b>1</b>	4	75	Undyed	36		CL-867

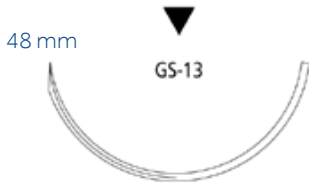

## REVERSE CUTTING

Needles Specifications	Size		Length cm	Color	Units Per box	Comments	Reorder code
	USP	EP					
<b>HOS-10</b> 1/2 Circle 26 mm 	<b>2-0</b>	3	3 × 45	Violet	12	D-Tach™	CL-356-MG
	<b>2-0</b>	3	3 × 45	Undyed	12	D-Tach™	CL-354-MG
	<b>0</b>	3.5	5 × 45	Undyed	12	D-Tach™	CL-55-MG
	<b>0</b>	3.5	3 × 45	Undyed	12	D-Tach™	CL-355-MG
	<b>1</b>	4	3 × 45	Violet	12	D-Tach™	CL-358-MG
	<b>1</b>	4	75	Undyed	36		CL-581
	<b>1</b>	4	75	Violet	36		CL-547
	<b>2</b>	5	75	Violet	36		CL-552
<b>C-25</b> 1/2 Circle 30 mm 	<b>2-0</b>	3	75	Violet	36		SL-622
	<b>0</b>	3.5	90	Undyed	36		SL-721
	<b>1</b>	4	90	Violet	36		SL-722
<b>GS-11</b> 1/2 Circle 37 mm 	<b>3-0</b>	2	75	Violet	36		CL-815
	<b>3-0</b>	2	75	Undyed	36		CL-542
	<b>2-0</b>	3	75	Violet	36		CL-816
	<b>2-0</b>	3	75	Undyed	36		CL-543
	<b>2-0</b>	3	90	Violet	36		CL-916
	<b>0</b>	3.5	75	Violet	36		CL-817
	<b>0</b>	3.5	75	Undyed	36		CL-544
	<b>0</b>	3.5	90	Violet	36		CL-917
	<b>0</b>	3.5	90	Undyed	36		CL-534
	<b>1</b>	4	75	Violet	36		CL-818
	<b>1</b>	4	75	Undyed	36		CL-545
	<b>1</b>	4	90	Violet	36		CL-918
	<b>1</b>	4	90	Undyed	36		CL-535
	<b>2</b>	5	90	Violet	36		CL-919
	<b>2</b>	5	90	Undyed	36		CL-533
	<b>GS-11 - GS-21</b> 1/2 Circle 37 mm Double-arm (Not to scale) 	<b>3-0</b>	2	150	Violet	36	
<b>2-0</b>		3	150	Violet	36		CL-992


## REVERSE CUTTING

Needles Specifications	Size		Length cm	Color	Units Per box	Comments	Reorder code
	USP	EP					
<b>HOS-11</b> 1/2 Circle 37 mm 	2-0	3	3 × 45	Undyed	12	D-Tach™	CL-374-MG
	2-0	3	3 × 45	Violet	12	D-Tach™	CL-369-MG
	2-0	3	75	Violet	36		CL-526
	2-0	3	75	Undyed	36		CL-536
	0	3.5	3 × 45	Undyed	12	D-Tach™	CL-375-MG
	0	3.5	3 × 45	Violet	12	D-Tach™	CL-370-MG
	0	3.5	75	Undyed	36		CL-537
	1	4	75	Undyed	36		CL-538
	1	4	75	Violet	36		CL-528
<b>GS-12</b> 1/2 Circle 40 mm 	2-0	3	75	Violet	36		CL-833
	2-0	3	75	Undyed	36		CL-595
	0	3.5	75	Violet	36		CL-804
	0	3.5	75	Undyed	36		CL-596
	0	3.5	90	Violet	36		CL-933
	1	4	75	Violet	36		CL-805
	1	4	75	Undyed	36		CL-597
	1	4	90	Violet	36		CL-936
	2	5	90	Violet	36		CL-937
<b>HOS-12</b> 1/2 Circle 40 mm 	2-0	3	3 × 45	Violet	12	D-Tach™	CL-382-MG
	0	3.5	3 × 45	Violet	12	D-Tach™	CL-380-MG
	0	3.5	3 × 45	Undyed	12	D-Tach™	CL-376-MG
	1	4	3 × 45	Violet	12	D-Tach™	CL-381-MG
	1	4	3 × 45	Undyed	12	D-Tach™	CL-377-MG
	1	4	5 × 75	Violet	12	D-Tach™	CL-97-MG
	2	5	3 × 45	Violet	12	D-Tach™	CL-383-MG




## REVERSE CUTTING

Needles Specifications	Size		Length cm	Color	Units Per box	Comments	Reorder code
	USP	EP					
<b>GS-13</b> 1/2 Circle    48 mm 	<b>0</b>	3.5	75	Violet	36		CL-568
	<b>0</b>	3.5	75	Undyed	36		CL-554
	<b>0</b>	3.5	90	Violet	36		CL-927
	<b>1</b>	4	75	Violet	36		CL-569
	<b>1</b>	4	75	Undyed	36		CL-555
	<b>1</b>	4	90	Violet	36		CL-928
	<b>2</b>	5	75	Violet	36		CL-570
	<b>2</b>	5	90	Undyed	36		CL-939
	<b>2</b>	5	90	Violet	36		CL-929
	<b>HOS-14</b> 1/2 Circle    55 mm 	<b>2</b>	5	75	Violet	36	



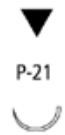


## PREMIUM REVERSE CUTTING

Needles Specifications	Size		Length cm	Color	Units Per box	Comments	Reorder code
	USP	EP					
<b>P-10</b> 3/8 Circle    11 mm 	<b>7-0</b>	0.5	45	Undyed	12		SL-1624G
	<b>6-0</b>	0.7	45	Undyed	36	Surgalloy™	SL-1625
	<b>6-0</b>	0.7	45	Undyed	12	Surgalloy™	SL1625G
	<b>5-0</b>	1	45	Undyed	36	Surgalloy™	SL-1689
	<b>5-0</b>	1	45	Undyed	12	Surgalloy™	SL-1689G




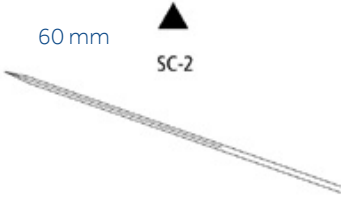
## PREMIUM REVERSE CUTTING

Needles Specifications	Size USP	EP	Length cm	Color	Units Per box	Comments	Reorder code
<b>P-13</b> 3/8 Circle 13 mm 	6-0	0.7	45	Undyed	36	Surgalloy™	SL-5688
	6-0	0.7	45	Undyed	12	Surgalloy™	SL-5688-G
	5-0	1	45	Undyed	12	Surgalloy™	SL-5687G
	5-0	1	45	Undyed	36	Surgalloy™	SL-5687
	4-0	1.5	45	Undyed	12	Surgalloy™	SL-5691G
	4-0	1.5	45	Undyed	36	Surgalloy™	SL-5690
<b>P-11</b> 3/8 Circle 16 mm 	5-0	1	45	Undyed	36	Surgalloy™	SL-5626
	4-0	1.5	45	Undyed	36	Surgalloy™	SL-5686
	4-0	1.5	45	Undyed	12	Surgalloy™	SL-5686G
	4-0	1.5	75	Undyed	36	Surgalloy™	SL-5620
	3-0	2	75	Undyed	36	Surgalloy™	SL-5621
<b>P-12</b> 3/8 Circle 19 mm 	5-0	1	45	Undyed	12	Surgalloy™	SL-5636-G
	5-0	1	45	Undyed	36	Surgalloy™	SL-5636
	4-0	1.5	45	Violet	36	Surgalloy™	SL-5629
	4-0	1.5	45	Undyed	12	Surgalloy™	SL-5627-G
	4-0	1.5	45	Undyed	36	Surgalloy™	SL-5627
	4-0	1.5	75	Undyed	36	Surgalloy™	SL-5637
	3-0	2	45	Undyed	12	Surgalloy™	SL-5628G
	3-0	2	45	Undyed	36	Surgalloy™	SL-5628
	3-0	2	75	Undyed	36	Surgalloy™	SL-5638


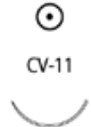




## PREMIUM REVERSE CUTTING

Needles Specifications	Size		Length cm	Color	Units Per box	Comments	Reorder code
	USP	EP					
<b>P-14</b> 3/8 Circle 24 mm 	4-0	1.5	45	Undyed	12	Surgalloy™	SL-5641G
	4-0	1.5	45	Undyed	36	Surgalloy™	SL-5641
	4-0	1.5	75	Undyed	36	Surgalloy™	SL-5678
	3-0	2	45	Undyed	12	Surgalloy™	SL-5640G
	3-0	2	45	Undyed	36	Surgalloy™	SL-5640
	3-0	2	75	Undyed	12	Surgalloy™	SL-5679G
	3-0	2	75	Undyed	36	Surgalloy™	SL-5679
	2-0	3	75	Undyed	36		SL-5680
<b>P-17</b> 3/8 Circle 30 mm 	2-0	3	75	Undyed	36		SL-1596
<b>P-21</b> 1/2 Circle 9 mm 	6-0	0.7	45	Undyed	12		SL-1510G
	5-0	1	45	Undyed	12		SL-1644G
	4-0	1.5	45	Undyed	12		SL-1645G
<b>P-22</b> 1/2 Circle 13 mm 	6-0	0.7	45	Undyed	36	Surgalloy™	SL-1613
	5-0	1	45	Undyed	36	Surgalloy™	SL-1614
	5-0	1	45	Undyed	12	Surgalloy™	SL-1614G
	4-0	1.5	45	Undyed	36	Surgalloy™	SL-1615
<b>P-24</b> 1/2 Circle 16 mm 	5-0	1	45	Undyed	12	Surgalloy™	SL-5632G
	5-0	1	45	Undyed	36	Surgalloy™	SL-5632
	4-0	1.5	45	Undyed	12	Surgalloy™	SL-5633G
	4-0	1.5	45	Undyed	36	Surgalloy™	SL-5633

## CONVENTIONAL CUTTING

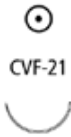
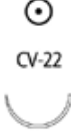



Needles Specifications	Size		Length cm	Color	Units Per box	Comments	Reorder code
	USP	EP					
<b>PC-13</b> 3/8 Circle 11 mm 	6-0	0.7	45	Undyed	36		SL-1660
	5-0	1	45	Undyed	36		SL-1659
	4-0	1.5	45	Undyed	36		SL-1658
<b>PC-11</b> 3/8 Circle 16 mm 	6-0	0.7	45	Undyed	12	Surgalloy™	SL-1655G
	5-0	1	45	Undyed	36		SL-1656
	4-0	1.5	45	Undyed	36		SL-1657
<b>PC-12</b> 3/8 Circle 19 mm 	4-0	1.5	45	Undyed	36		SL-1653
	3-0	2	45	Undyed	36		SL-1654
<b>SC-2</b> Straight 60 mm 	4-0	1.5	75	Violet	36		SL-654
	4-0	1.5	75	Undyed	36		SL-644
	3-0	2	75	Violet	36		SL-653
	3-0	2	75	Undyed	36		SL-643
	2-0	3	75	Violet	36		SL-652
	2-0	3	75	Undyed	36		SL-642

## TAPER POINT





Needles Specifications	Size USP	EP	Length cm	Color	Units Per box	Comments	Reorder code
<b>MV-135-5</b> 3/8 Circle 6 mm 	8-0	0.4	13	Violet	12	Surgalloy™	L-2800
<b>CV-11</b> 3/8 Circle 13 mm 	6-0	0.7	75	Violet	36		GL-889-2
	5-0	1	75	Violet	36		GL-890
<b>CV-15</b> 3/8 Circle 17 mm 	4-0	1.5	75	Violet	36		GL-881
	3-0	2	75	Violet	36		GL-882
<b>CV-13</b> 3/8 Circle 22 mm 	4-0	1.5	75	Violet	36		CL-304
	3-0	2	75	Violet	36		GL-303
	2-0	3	75	Violet	36		GL-302
<b>CV-18</b> 3/8 Circle 26 mm 	2-0	3	75	Violet	36		CL-457
<b>CV-17</b> 3/8 Circle 32 mm 	3-0	2	75	Violet	36		GL-522
	2-0	3	75	Violet	36		GL-523
	0	3.5	75	Violet	36		GL-524





## TAPER POINT

Needles Specifications	Size USP	EP	Length cm	Color	Units Per box	Comments	Reorder code
<b>CVF-21</b> 1/2 Circle 12 mm 	6-0	0.7	45	Violet	36	Surgalloy™	UL-102
	6-0	0.7	45	Undyed	36	Surgalloy™	UL-101
	5-0	1	45	Undyed	36	Surgalloy™	SL-432
	4-0	1.5	45	Undyed	36		SL-431
<b>CV-22</b> 1/2 Circle 13 mm 	5-0	1	5 × 45	Violet	12	D-Tach™	GL-45-MG
	5-0	1	75	Undyed	36		GL-885
	4-0	1.5	5 × 45	Violet	12	D-Tach™	GL-44-MG
	4-0	1.5	75	Undyed	36		GL-884
	4-0	1.5	75	Violet	36		CL-782
<b>CV-23</b> 1/2 Circle 17 mm 	6-0	0.7	75	Violet	36		UL-201-2
	5-0	1	75	Violet	36		UL-202
	5-0	1	75	Undyed	36		UL-213
	4-0	1.5	5 × 45	Violet	12	D-Tach™	GL-46-MG
	4-0	1.5	75	Violet	36	Surgalloy™	UL-303
	4-0	1.5	75	Undyed	36		UL-212
	3-0	2	5 × 45	Undyed	12	D-Tach™	UL-47-MG
	3-0	2	5 × 45	Violet	12	D-Tach™	GL-47-MG
	3-0	2	75	Violet	36		UL-204
	3-0	2	75	Undyed	36		UL-215
	2-0	3	75	Undyed	36		UL-216
	2-0	3	75	Violet	36		UL-205
<b>CVF-23</b> 1/2 Circle 17 mm 	4-0	1.5	75	Undyed	36	Surgalloy™	UL-214
	4-0	1.5	75	Violet	36	Surgalloy™	UL-203
<b>CVF-23</b> 1/2 Circle 17 mm Double-arm 	4-0	1.5	90	Violet	36	Surgalloy™	CL-219


## TAPER POINT

Needles Specifications	Size USP	EP	Length cm	Color	Units Per box	Comments	Reorder code
<b>CV-24</b> 1/2 Circle 20 mm 	4-0	1.5	75	Violet	36		CL-460
	3-0	2	75	Violet	36		CL-461
	2-0	3	75	Violet	36		CL-462
<b>CV-25</b> 1/2 Circle 22 mm 	4-0	1.5	5 × 45	Violet	12	D-Tach™	GL-34-MG
	4-0	1.5	75	Violet	36		GL-181
	4-0	1.5	75	Undyed	36		GL-191
	3-0	2	5 × 45	Violet	12	D-Tach™	GL-33-MG
	3-0	2	75	Violet	36		GL-182
	3-0	2	75	Undyed	36		GL-192
	3-0	2	5 × 75	Violet	12	D-Tach™	GL-182-M
	3-0	2	120	Violet	24	Surgalloy™ Endoscopy	EL-440
	2-0	3	5 × 45	Violet	12	D-Tach™	GL-32-MG
	2-0	3	75	Violet	36		GL-183
<b>CV-25</b> 1/2 Circle 22 mm Double-arm 	3-0	2	90	Violet	36		GL-282
<b>GS-23</b> 1/2 Circle 22 mm 	3-0	2	75	Violet	36		CL-892
	2-0	3	5 × 45	Undyed	12	D-Tach™	CL-42-MG
	2-0	3	75	Violet	36		CL-893
	0	3.5	75	Violet	36		CL-894

## TAPER POINT

Needles Specifications	Size		Length cm	Color	Units Per box	Comments	Reorder code	
	USP	EP						
<b>V-20</b> 1/2 Circle 26 mm 	5-0	1	75	Violet	36		GL-120	
	4-0	1.5	5 × 45	Violet	12	D-Tach™	GL-64-MG	
	4-0	1.5	75	Undyed	36	Surgalloy™	GLS-321	
	4-0	1.5	75	Violet	36	Surgalloy™	GLS-121	
	4-0	1.5	5 × 75	Violet	12	D-Tach™	GL-68-M	
	4-0	1.5	75	Undyed	36		GL-321	
	4-0	1.5	75	Violet	36		GL-121	
	3-0	2	5 × 45	Undyed	12	D-Tach™	GL-69-MG	
	3-0	2	5 × 45	Violet	12	D-Tach™	GL-63-MG	
	3-0	2	75	Undyed	36	Surgalloy™	GLS-322	
	3-0	2	75	Violet	36	Surgalloy™	GLS-122	
	3-0	2	5 × 75	Violet	12	D-Tach™	GL-67-M	
	3-0	2	75	Undyed	36		GL-322	
	3-0	2	75	Violet	36		GL-122	
	3-0	2	90	Violet	36		GL-226	
	2-0	3	5 × 45	Violet	12	D-Tach™	GL-62-MG	
	2-0	3	5 × 45	Undyed	12	D-Tach™	GL-61-MG	
	2-0	3	75	Undyed	36	Surgalloy™	GLS-323	
	2-0	3	75	Violet	36	Surgalloy™	GLS-123	
	2-0	3	5 × 75	Violet	12	D-Tach™	GL-66-M	
	2-0	3	75	Undyed	36		GL-323	
	2-0	3	75	Violet	36		GL-123	
	2-0	3	90	Violet	36		GL-125	
	0	3.5	75	Undyed	36		GL-324	
	0	3.5	75	Violet	36		GL-124	
	<b>V-20</b> 1/2 Circle 26 mm Double-arm 	3-0	2	90	Violet	36		GL-332
		2-0	3	90	Violet	36		GL-333

## TAPER POINT

Needles Specifications	Size		Length cm	Color	Units Per box	Comments	Reorder code
	USP	EP					
<b>GS-22</b> 1/2 Circle 27 mm 	3-0	2	75	Violet	36		CL-882
	3-0	2	75	Undyed	36		CL-829
	2-0	3	5 × 45	Violet	12	D-Tach™	CL-16-MG
	2-0	3	75	Undyed	36		CL-830
	2-0	3	75	Violet	36		CL-883
	0	3.5	5 × 45	Violet	12	D-Tach™	CL-17-MG
	0	3.5	75	Violet	36	Surgalloy™	CLS-884
	0	3.5	5 × 75	Violet	12	D-Tach™	CL-886-M
	0	3.5	75	Violet	36		CL-884
	0	3.5	75	Undyed	36		CL-831
	0	3.5	90	Violet	36		CL-330
	1	4	75	Violet	36		CL-885
	1	4	5 × 75	Violet	12	D-Tach™	CL-885-M
	1	4	75	Undyed	36		CL-832
	1	4	90	Violet	36		CL-331
	2	5	75	Violet	36		CL-232
	<b>V-30</b> 1/2 Circle 30 mm 	3-0	2	75	Violet	36	
2-0		3	75	Violet	36		GL-223
2-0		3	75	Undyed	36		GL-228
2-0		3	5 × 75	Violet	12	D-Tach™	GL-223-M
0		3.5	75	Violet	36		GL-224
1		4	75	Violet	36		GL-225

## TAPER POINT

Needles Specifications	Size		Length cm	Color	Units Per box	Comments	Reorder code
	USP	EP					
<b>GS-21</b> 1/2 Circle 37 mm 	4-0	1.5	75	Violet	36		CL-809
	3-0	2	5 × 45	Violet	12	D-Tach™	CL-8-MG
	3-0	2	75	Undyed	36		CL-842
	3-0	2	75	Violet	36		CL-810
	3-0	2	90	Violet	36		CL-922
	3-0	2	90	Undyed	36		CL-952
	2-0	3	5 × 45	Violet	12	D-Tach™	CL-9-MG
	2-0	3	5 × 45	Undyed	12	D-Tach™	CL-13-MG
	2-0	3	75	Undyed	36		CL-843
	2-0	3	75	Violet	36		CL-811
	2-0	3	5 × 75	Violet	12	D-Tach™	CL-4M
	2-0	3	90	Violet	36		CL-923
	2-0	3	90	Undyed	36		CL-953
	0	3.5	5 × 45	Undyed	12	D-Tach™	CL-14-MG
	0	3.5	5 × 45	Violet	12	D-Tach™	CL-10-MG
	0	3.5	75	Undyed	36		CL-844
	0	3.5	75	Violet	36		CL-812
	0	3.5	5 × 75	Violet	12	D-Tach™	CL-5-M
	0	3.5	5 × 75	Undyed	12	D-Tach™	CL-2M
	0	3.5	75	Violet	24		CL-812
	0	3.5	90	Undyed	36		CL-954
	0	3.5	90	Violet	36		CL-924
	1	4	5 × 45	Undyed	12	D-Tach™	CL-15-MG
	1	4	5 × 45	Violet	12	D-Tach™	CL-11-MG
	1	4	75	Undyed	36		CL-845
	1	4	75	Violet	36		CL-813
	1	4	5 × 75	Violet	12	D-Tach™	CL-6-M
	1	4	90	Violet	36		CL-925
	1	4	90	Undyed	36		CL-955
	2	5	75	Violet	36		CL-814
2	5	90	Violet	36		CL-926	

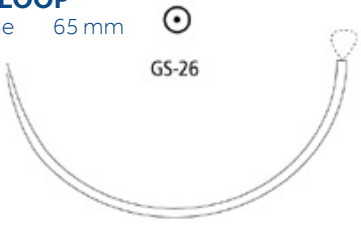




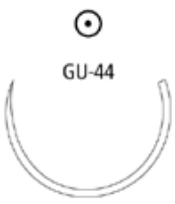
## TAPER POINT

Needles Specifications	Size		Length cm	Color	Units Per box	Comments	Reorder code
	USP	EP					
<b>V-26</b> 1/2 Circle 37 mm 	<b>3-0</b>	2	75	Violet	36		GL-126
	<b>2-0</b>	3	75	Violet	36		GL-127
	<b>0</b>	3.5	75	Violet	36		GL-128
	<b>1</b>	4	75	Violet	36		GL-129
<b>GS-24</b> 1/2 Circle 40 mm 	<b>3-0</b>	2	75	Violet	36		CL-800
	<b>3-0</b>	2	90	Undyed	36		CL-962
	<b>3-0</b>	2	90	Violet	36		CL-912
	<b>2-0</b>	3	5 × 45	Violet	12	D-Tach™	CL-22-MG
	<b>2-0</b>	3	75	Undyed	36		CL-853
	<b>2-0</b>	3	75	Violet	36		CL-801
	<b>2-0</b>	3	90	Violet	36		CL-913
	<b>2-0</b>	3	90	Undyed	36		CL-963
	<b>0</b>	3.5	5 × 45	Violet	12	D-Tach™	CL-20-MG
	<b>0</b>	3.5	75	Violet	36		CL-802
	<b>0</b>	3.5	90	Violet	36		CL-914
	<b>0</b>	3.5	90	Undyed	36		CL-964
	<b>1</b>	4	5 × 45	Violet	12	D-Tach™	CL-21-MG
	<b>1</b>	4	75	Violet	36		CL-803
	<b>1</b>	4	75	Undyed	36		CL-855
	<b>1</b>	4	90	Violet	36		CL-915
	<b>1</b>	4	90	Undyed	36		CL-965
	<b>2</b>	5	5 × 45	Violet	6	D-Tach™	CL-623-M
	<b>2</b>	5	5 × 60	Violet	12	D-Tach™	CL-622-M
<b>2</b>	5	75	Violet	36		CL-779	
<b>2</b>	5	90	Violet	36		CL-911	
<b>GS-24</b> 1/2 Circle 40 mm Double-arm 	<b>3-0</b>	2	150	Violet	36		CL-990

## TAPER POINT



Needles Specifications	Size		Length cm	Color	Units Per box	Comments	Reorder code
	USP	EP					
<b>GS-25</b> 1/2 Circle 48 mm 	<b>3-0</b>	2	75	Violet	36		CL-862
	<b>3-0</b>	2	90	Undyed	36		CL-972
	<b>2-0</b>	3	5 × 45	Undyed	12	D-Tach™	CL-36-MG
	<b>2-0</b>	3	75	Violet	36		CL-863
	<b>2-0</b>	3	90	Undyed	36		CL-973
	<b>2-0</b>	3	90	Violet	36		CL-903
	<b>0</b>	3.5	5 × 45	Undyed	12	D-Tach™	CL-35-MG
	<b>0</b>	3.5	5 × 45	Violet	12	D-Tach™	CL-30-MG
	<b>0</b>	3.5	75	Violet	36		CL-864
	<b>0</b>	3.5	90	Undyed	36		CL-974
	<b>0</b>	3.5	90	Violet	36		CL-904
	<b>1</b>	4	5 × 45	Violet	12	D-Tach™	CL-31-MG
	<b>1</b>	4	75	Violet	36		CL-865
	<b>1</b>	4	75	Violet	12	D-Tach™	CL-64-M
	<b>1</b>	4	90	Undyed	36		CL-975
	<b>1</b>	4	90	Violet	36		CL-905
	<b>2</b>	5	5 × 45	Violet	6	D-Tach™	CL-633-M
	<b>2</b>	5	5 × 60	Violet	12	D-Tach™	CL-644-M
	<b>2</b>	5	75	Violet	36		CL-866
	<b>2</b>	5	5 × 75	Violet	6	D-Tach™	CL-665-M
<b>2</b>	5	90	Violet	36		CL-906	
<b>GS-26</b> 1/2 Circle 65 mm 	<b>2-0</b>	3	75	Undyed	24		CL-63
	<b>1</b>	4	5 × 75	Violet	12		CL-60-M
	<b>1</b>	4	90	Violet	36		CL-968
	<b>1</b>	4	150	Violet	24		CL-824
	<b>2</b>	5	5 × 75	Violet	6		CL-659-M
	<b>2</b>	5	100	Undyed	36		CL-969
	<b>2</b>	5	150	Undyed	24		CL-880
	<b>2</b>	5	150	Violet	24		CL-825

## TAPER POINT




Needles Specifications	Size		Length cm	Color	Units Per box	Comments	Reorder code
	USP	EP					
<b>GS-26 LOOP</b> 1/2 Circle 65 mm 	<b>1</b>	4	120	Violet	24	Loop	CL-824L
	<b>2</b>	5	240	Violet	24	Loop	CL-835-L
<b>GS-27</b> 1/2 Circle 76 mm 	<b>2-0</b>	3	75	Undyed	24		CL-53
	<b>1</b>	4	75	Violet	24	455	SL-603
<b>GU-46</b> 5/8 Circle 27 mm 	<b>3-0</b>	2	75	Violet	36		UL-879
	<b>2-0</b>	3	75	Violet	36		UL-878
	<b>0</b>	3.5	75	Violet	36		UL-877
	<b>1</b>	4	75	Violet	36		UL-876
<b>GU-46</b> 5/8 Circle 27 mm Double-arm 	<b>0</b>	3.5	90	Violet	36		UL-240
<b>GU-45</b> 5/8 Circle 37 mm 	<b>2-0</b>	3	75	Violet	36		UL-245
	<b>0</b>	3.5	75	Violet	36		UL-246
<b>GU-44</b> 5/8 Circle 40 mm 	<b>0</b>	3.5	75	Violet	36		UL-146





## TAPER POINT

Needles Specifications			Size		Length cm	Color	Units Per box	Comments	Reorder code
			USP	EP					
<b>SK</b> SKI	23 mm		3-0	2	18	Violet	36	Endoscopy	EL415
			3-0	2	120	Violet	24	Endoscopy	EL450
<b>TJ-32</b> 1/2 Circle	32 mm		1	4	75	Violet	12		CL-673G
			1	4	75	Violet	36		CL-673



## HEAVY TAPER POINT

Needles Specifications			Size		Length cm	Color	Units Per box	Comments	Reorder code
			USP	EP					
<b>HGS-22</b> 1/2 Circle	27 mm		2-0	3	5 × 45	Violet	12	D-Tach™	CL-86-MG
			0	3.5	5 × 45	Undyed	12	D-Tach™	CL-87-MG
			1	4	75	Violet	36		CL-871
<b>HGS-20</b> 1/2 Circle	34 mm		2-0	3	5 × 75	Violet	12	D-Tach™	CL-996-MG
			0	3.5	5 × 75	Violet	12	D-Tach™	CL-997-MG
			1	4	5 × 75	Violet	12	D-Tach™	CL-998-MG
<b>HGS-21</b> 1/2 Circle	37 mm		2-0	3	5 × 45	Violet	12	D-Tach™	CL-202-M
			2-0	3	75	Violet	36		CL-930
			2-0	3	90	Violet	36		CL-970
			0	3.5	5 × 45	Violet	12	D-Tach™	CL-200-M
			0	3.5	75	Violet	36		CL-931
			0	3.5	5 × 75	Violet	12	D-Tach™	CL-71-M
			0	3.5	90	Violet	36		CL-971
			1	4	5 × 45	Violet	12	D-Tach™	CL-201-M
			1	4	75	Violet	36		CL-932
			1	4	90	Violet	36		CL-239
			2	5	75	Violet	36		CL-934




## HEAVY TAPER POINT

Needles Specifications	Size		Length cm	Color	Units Per box	Comments	Reorder code
	USP	EP					
<b>HGS-21 LOOP</b> 1/2 Circle 37 mm 	0	3.5	150	Violet	24	Loop	CL-871-L
<b>HGU-46</b> 5/8 Circle 27 mm 	0	3.5	45	Violet	36		CL-576
	0	3.5	75	Violet	36		CL-677
	1	4	75	Violet	36		UL-676



## TAPERCUTTING

Needles Specifications	Size		Length cm	Color	Units Per box	Comments	Reorder code
	USP	EP					
<b>KV-7</b> 1/2 Circle 26 mm 	3-0	2	75	Violet	36	Surgalloy™	CL-510
	2-0	3	75	Violet	36	Surgalloy™	CL-511
<b>KV-34</b> 1/2 Circle 37 mm 	3-0	2	90	Violet	36		CL-942
	2-0	3	75	Violet	36		CL-837
	2-0	3	90	Violet	36		CL-943
	2-0	3	90	Undyed	36		CL-946
	0	3.5	75	Violet	36		CL-838
	0	3.5	90	Undyed	36		CL-515
	0	3.5	90	Violet	36		CL-944
	1	4	75	Violet	36		CL-839
	1	4	90	Undyed	36		CL-516
	1	4	90	Violet	36		CL-945
	2	5	90	Violet	36		CL-517




## TAPERCUTTING

Needles Specifications	Size		Length cm	Color	Units Per box	Comments	Reorder code
	USP	EP					
<b>KV-37</b> 1/2 Circle 40 mm 	2-0	3	75	Violet	36		CL-603
	0	3.5	75	Violet	36		CL-602
<b>KV-40</b> 1/2 Circle 48 mm 	0	3.5	90	Violet	36		CL-739
	1	4	75	Violet	36		CL-840
	1	4	90	Violet	36		CL-940
	2	5	75	Violet	36		CL-841
	2	5	90	Violet	36		CL-941
<b>KV-40 LOOP</b> 1/2 Circle 48 mm 	1	4	150	Violet	24	Loop	CL-740-L
	2	5	150	Violet	24	Loop	CL-741-L


## BLUNT POINT

Needles Specifications	Size		Length cm	Color	Units Per box	Comments	Reorder code
	USP	EP					
<b>BTV-20</b> 1/2 Circle 26 mm 	3-0	2	5 × 45	Undyed	12	D-Tach™, ProtectPoint™	GLT-69-MG
<b>BGS-21</b> 1/2 Circle 37 mm 	2-0	3	90	Violet	36	ProtectPoint™	CL-956
	2-0	3	90	Undyed	36	ProtectPoint™	CL-259
	0	3.5	5 × 45	Violet	12	D-Tach™, ProtectPoint™	CL-250-M
	0	3.5	5 × 45	Undyed	12	D-Tach™, ProtectPoint™	CL-240-M
	0	3.5	90	Violet	36	ProtectPoint™	CL-957
	1	4	90	Violet	36	ProtectPoint™	CL-958

## BLUNT POINT

Needles Specifications	Size		Length cm	Color	Units Per box	Comments	Reorder code
	USP	EP					
<b>BTP-1</b> 1/2 Circle 37 mm 	<b>3-0</b>	2	90	Undyed	36	ProtectPoint™	CLT-952
	<b>2-0</b>	3	5 × 45	Violet	12	D-Tach™, ProtectPoint™	CLT-9-MG
	<b>2-0</b>	3	5 × 45	Undyed	12	D-Tach™, ProtectPoint™	CLT-13-MG
	<b>2-0</b>	3	75	Undyed	36	ProtectPoint™	CLT-259
	<b>2-0</b>	3	90	Undyed	36	ProtectPoint™	CLT-953
	<b>0</b>	3.5	5 × 45	Violet	12	D-Tach™, ProtectPoint™	CLT-250-MG
	<b>0</b>	3.5	75	Violet	36		CLT-957
	<b>0</b>	3.5	90	Undyed	36	ProtectPoint™	CLT-954
	<b>1</b>	4	5 × 45	Undyed	12	D-Tach™, ProtectPoint™	CLT-15-MG
	<b>1</b>	4	5 × 45	Violet	12	D-Tach™, ProtectPoint™	CLT-11-MG
	<b>1</b>	4	75	Violet	36	ProtectPoint™	CLT-958
	<b>1</b>	4	90	Undyed	36	ProtectPoint™	CLT-955
	<b>BGS-24</b> 1/2 Circle 40 mm 	<b>0</b>	3.5	90	Violet	36	ProtectPoint™
<b>1</b>		4	90	Violet	36	ProtectPoint™	CL-951
<b>BTP-X</b> 1/2 Circle 48 mm 	<b>0</b>	3.5	5 × 45	Violet	12	D-Tach™, ProtectPoint™	CLT-850-MG
	<b>0</b>	3.5	5 × 45	Undyed	12	D-Tach™, ProtectPoint™	CLT-20-MG
	<b>1</b>	4	5 × 45	Violet	12	D-Tach™, ProtectPoint™	CLT-851-MG
	<b>1</b>	4	5 × 45	Undyed	12	D-Tach™, ProtectPoint™	CLT-21-MG
	<b>1</b>	4	90	Undyed	36	ProtectPoint™	CLT-855
	<b>2</b>	5	5 × 45	Violet	6	D-Tach™, ProtectPoint™	CLT-623-M




## BLUNT POINT

Needles Specifications	Size		Length cm	Color	Units Per box	Comments	Reorder code
	USP	EP					
<b>BGS-25</b> 1/2 Circle    48 mm 	<b>0</b>	3.5	90	Violet	36	ProtectPoint™	CL-948
	<b>1</b>	4	90	Violet	36	ProtectPoint™	CL-949
<b>BGS-29</b> 1/2 Circle    63 mm 	<b>2</b>	5	2 × 90	Violet	12	ProtectPoint™	2-CL-852
<b>BP-27</b> 1/2 Circle    64 mm 	<b>1</b>	4	90	Violet	36	ProtectPoint™	CL-271
	<b>2</b>	5	90	Violet	36	ProtectPoint™	CL-270
<b>BGST-29</b> 1/2 Circle    65 mm 	<b>2</b>	5	5 × 75	Violet	6	D-Tach™, ProtectPoint™	CLT-633-M

## DIAMOND POINT




Needles Specifications	Size		Length cm	Color	Units Per box	Comments	Reorder code
	USP	EP					
<b>KV-25</b> 1/2 Circle    22 mm 	<b>1</b>	4	75	Undyed	36		CL-513
<b>KV-30</b> 1/2 Circle    27 mm 	<b>0</b>	3.5	90	Undyed	36		CL-730

## MULTIPURPOSE



Needles Specifications	Size		Length cm	Color	Units Per box	Comments	Reorder code
	USP	EP					
<b>3.16</b> 3/8 Circle    16 mm 	<b>4-0</b>	1.5	45	Undyed	36		CL-005
<b>4.21 (H)*</b> 1/2 Circle    20 mm 	<b>2-0</b>	3	75	Violet	36		CL-015
<b>4.31 (H)*</b> 1/2 Circle    30 mm 	<b>0</b>	3.5	75	Violet	36		CL-030

**\*(H) = Heavy**

## MULTIPURPOSE

Needles Specifications	Size		Length cm	Color	Units Per box	Comments	Reorder code
	USP	EP					
<b>4.41</b> 1/2 Circle 40 mm 	<b>2-0</b>	3	75	Violet	36		CL-045
	<b>0</b>	3.5	75	Violet	36		UL-046
	<b>1</b>		75	Violet	36		CL-047
<b>4.5</b> 1/2 Circle 50 mm (Not to scale) 	<b>2</b>	5	90	Violet	36		CL-055
<b>4.51 (H)</b> 1/2 Circle 50 mm 	<b>2-0</b>	3	75	Violet	36		CL-060

## SPATULA





Needles Specifications	Size		Length cm	Color	Units Per box	Comments	Reorder code
	USP	EP					
<b>SE-100-8</b> 1/4 Circle 6.3 mm Double-arm Curvature 100° 	<b>6-0</b>	0.7	45	Violet	12	Surgalloy™	L-1750K
<b>SS-29</b> 1/4 Circle 7.4 mm Double-arm Curvature 110° 	<b>6-0</b>	0.7	30	Violet	12	Surgalloy™	L-1756K
	<b>6-0</b>	0.7	45	Violet	12	Surgalloy™	L-1755K

## SPATULA

Needles Specifications	Size USP	EP	Length cm	Color	Units Per box	Comments	Reorder code	
<b>SS-14</b> 1/4 Circle 8.6 mm Double-arm Curvature 115°		6-0	0.7	45	Violet	12	Surgalloy™	L-1791K
		6-0	0.7	45	Undyed	12	Surgalloy™	L-1670-K
		5-0	1	45	Violet	12	Surgalloy™	L-1792K
		5-0	1	45	Undyed	12	Surgalloy™	L-1669-K
<b>SS-24</b> 1/4 Circle 8.6 mm Double-arm Curvature 95°		6-0	0.7	20	Violet	12		L-1740-K
		6-0	0.7	30	Violet	12	Surgalloy™	L-1742K
		6-0	0.7	30	Undyed	12	Surgalloy™	L-1745K
		6-0	0.7	45	Undyed	24	Surgalloy™	L-2752K
		5-0	1	20	Violet	12	Surgalloy™	L-1739-K
		5-0	1	30	Violet	12	Surgalloy™	L-1785K
		5-0	1	45	Undyed	24	Surgalloy™	L-2753K
<b>SE-140-8</b> 3/8 Circle 6.3 mm Curvature 140°		8-0	0.4	30	Violet	12	Surgalloy™	L-2750K
		7-0	0.5	30	Violet	12	Surgalloy™	L-2751K
<b>SE-140-8</b> 3/8 Circle 6.3 mm Double-arm Curvature 140°		8-0	0.4	20	Violet	12	Surgalloy™	L-2747K
		8-0	0.4	30	Violet	12	Surgalloy™	L-2748K
		8-0	0.4	45	Violet	12	Surgalloy™	L-5118K
		7-0	0.5	30	Violet	12	Surgalloy™	L-1796K
		7-0	0.5	45	Violet	12	Surgalloy™	L-1897K
<b>SS-28</b> 3/8 Circle 7.4 mm Double-arm Curvature 145°		6-0	0.7	45	Violet	12	Surgalloy™	L-1562-K



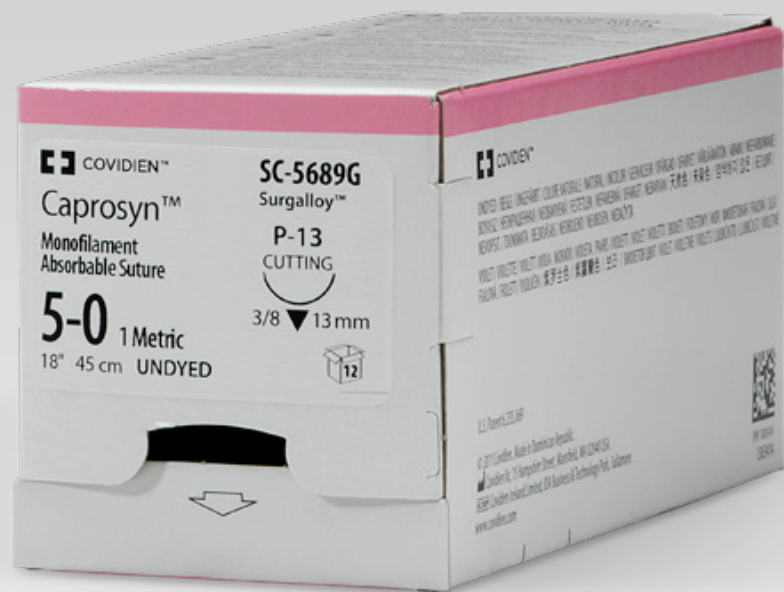
## SPATULA

Needles Specifications	Size USP	EP	Length cm	Color	Units Per box	Comments	Reorder code
<b>SE-160-8</b> 1/2 Circle 5.5 mm Double-arm Curvature 160°							
							
	<b>SE-160-8</b>						
							
<b>SS-2</b> 1/2 Circle 8.6 mm Double-arm Curvature 200°							
							
	<b>SS-2</b>						
							

# CAPROSYN™

## Monofilament Absorbable Suture

There when you need it. Gone when you don't.







Please refer to the individual package IFUs for complete instructions, indications, contraindications, warnings and precautions.

<b>Description</b>	Caprosyn™ monofilament synthetic absorbable sutures are prepared from Polyglytone™ 6211 synthetic polyester which is composed of glycolide, caprolactone, trimethylene carbonate, and lactide.
<b>Indications</b>	Caprosyn™ synthetic absorbable sutures are indicated for use in general soft tissue approximation and/or ligation, but not for use in cardiovascular or neurological surgery, microsurgery, or ophthalmic surgery.
<b>Actions</b>	Caprosyn™ synthetic absorbable sutures elicit a minimal acute inflammatory reaction in tissues, which is followed by gradual encapsulation of the suture by fibrous connective tissue. Progressive loss of tensile strength and eventual absorption of Caprosyn™ synthetic absorbable sutures occurs by means of hydrolysis.
<b>Tensile Strength</b>	Implantation studies in animals indicate that Caprosyn™ suture retains a minimum of 50 to 60% USP and EP knot strength at 5 days post implantation, and a minimum of 20 – 30% of USP and EP knot strength at 10 days post implantation. All of the original tensile strength is essentially lost by 21 days post implantation.
<b>Sterilization</b>	Ethylene oxide.
<b>Absorbtion</b>	Absorption begins as a loss of tensile strength followed by a loss of mass. Absorption of Caprosyn™ synthetic absorbable suture is essentially complete by 56 days.
<b>Packaging</b>	Caprosyn™ synthetic absorbable sutures are available undyed (natural) or dyed with D&C Violet No. 2 in USP sizes 1 (4 metric) through 6–0 (0.7 Metric). The sutures are supplied sterile, in pre-cut lengths or ligating reels, non-needled or affixed to various needle types. The sutures are available in box quantities of one, two and three dozen.


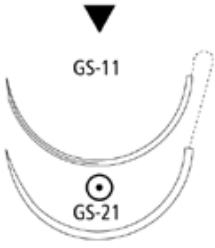
# Caprosyn™

Caprosyn™ suture is a brand you have trusted for over 10 years with a proven history of performance.



REEL							
Needles Specifications	Size		Length cm	Color	Units Per box	Comments	Reorder code
	USP	EP					
	3-0	2	150	Undyed	24		C-212
	0	3.5	150	Undyed	36		C-214

REVERSE CUTTING							
Needles Specifications	Size		Length cm	Color	Units Per box	Comments	Reorder code
	USP	EP					
<b>C-13</b> 3/8 Circle 19 mm 	5-0	1	45	Undyed	36		SC-3634
	4-0	1.5	45	Undyed	12		SC-690G
	4-0	1.5	75	Undyed	36		SC-691
	3-0	2	45	Undyed	12		SC-694G
	3-0	2	75	Undyed	36		SC-692
	2-0	3	75	Undyed	36		SC-693
<b>C-14</b> 3/8 Circle 24 mm 	3-0	2	75	Undyed	12	Surgalloy™	SC-693G
<b>C-15</b> 3/8 Circle 26 mm 	3-0	2	75	Undyed	36		SC-587
	2-0	3	75	Undyed	36		SC-586
<b>C-16</b> 3/8 Circle 30 mm 	0	3.5	75	Undyed	36	Surgalloy™	SC-585





## REVERSE CUTTING

Needles Specifications	Size		Length cm	Color	Units Per box	Comments	Reorder code
	USP	EP					
<b>GS-11</b> 1/2 Circle 37 mm 	2-0	3	75	Undyed	36		CC-543
<b>GS-11</b> 1/2 Circle 37 mm  <b>GS-21</b> Taper Point 1/2 Circle 37 mm Double-arm 	2-0	3	150	Undyed	36		CC-983





## PREMIUM REVERSE CUTTING

Needles Specifications	Size		Length cm	Color	Units Per box	Comments	Reorder code
	USP	EP					
<b>P-10</b> 3/8 Circle 11 mm 	6-0	0.7	45	Undyed	12	Surgalloy™	SC-3625G
	5-0	1	45	Undyed	12	Surgalloy™	SC-3626G
<b>P-13</b> 3/8 Circle 13 mm 	6-0	0.7	45	Undyed	12	Surgalloy™	SC-5688G
	6-0	0.7	45	Violet	12	Surgalloy™	SC-5588G
	5-0	1	45	Undyed	12	Surgalloy™	SC-5689G
	5-0	1	45	Violet	12	Surgalloy™	SC-5587G
	4-0	1.5	45	Undyed	12	Surgalloy™	SC-5690G
	4-0	1.5	45	Violet	12	Surgalloy™	SC-5590G





## PREMIUM REVERSE CUTTING

Needles Specifications	Size USP	EP	Length cm	Color	Units Per box	Comments	Reorder code
<b>P-11</b> 3/8 Circle 16 mm 	6-0	0.7	45	Undyed	12	Surgalloy™	SC-5635G
	5-0	1	45	Undyed	12	Surgalloy™	SC-5626G
	5-0	1	45	Undyed	36	Surgalloy™	SC-5626
	4-0	1.5	45	Undyed	12	Surgalloy™	SC-5686G
<b>P-12</b> 3/8 Circle 19 mm 	5-0	1	45	Undyed	12	Surgalloy™	SC-5616G
	4-0	1.5	45	Undyed	12	Surgalloy™	SC-5617G
	4-0	1.5	45	Violet	12	Surgalloy™	SC-5627G
	4-0	1.5	75	Undyed	12	Surgalloy™	SC-5637-G
	3-0	2	45	Undyed	12	Surgalloy™	SC-5618G
	3-0	2	75	Undyed	12	Surgalloy™	SC-5638-G
<b>P-14</b> 3/8 Circle 24 mm 	4-0	1.5	45	Undyed	36	Surgalloy™	SC-5679
	4-0	1.5	75	Undyed	36	Surgalloy™	SC-5779
	3-0	2	45	Undyed	36	Surgalloy™	SC-5640
	3-0	2	75	Undyed	12	Surgalloy™	SC-5780G
<b>P-24</b> 1/2 Circle 16 mm 	4-0	1.5	45	Undyed	12	Surgalloy™	SC-5643G


## CONVENTIONAL CUTTING

Needles Specifications	Size		Length cm	Color	Units Per box	Comments	Reorder code
	USP	EP					
<b>SC-1</b> Straight    51 mm  	3-0	2	75	Undyed	36	Surgalloy™	SC-543
	2-0	3	75	Undyed	36	Surgalloy™	SC-542
<b>SC-2</b> Straight    60 mm  	4-0	1.5	75	Undyed	36	Surgalloy™	SC-644
	3-0	2	75	Undyed	36		SC-643


## DERMAX PREMIUM X CUTTING

Needles Specifications	Size		Length cm	Color	Units Per box	Comments	Reorder code
	USP	EP					
<b>DX-13</b> 3/8 Circle    13 mm  DX-13 	5-0	1	45	Undyed	12	DermaX	XC-834
<b>DX-16</b> 3/8 Circle    16 mm  DX-16 	5-0	1	45	Undyed	12	DermaX	XC-844
	4-0	1.5	45	Undyed	12	DermaX	XC-845

## DERMAX PREMIUM X CUTTING

Needles Specifications	Size		Length cm	Color	Units Per box	Comments	Reorder code
	USP	EP					
<b>DX-19</b> 3/8 Circle    19 mm 	4-0	1.5	45	Undyed	12	DermaX	XC-823

## TAPER POINT






Needles Specifications	Size		Length cm	Color	Units Per box	Comments	Reorder code
	USP	EP					
<b>CV-23</b> 1/2 Circle    17 mm 	5-0	1	75	Undyed	36	Surgalloy™	UC-202
	4-0	1.5	75	Undyed	36	Surgalloy™	UC-203
	4-0	1.5	75	Violet	36	Surgalloy™	UC-214
	3-0	2	75	Violet	36		UC-213
	3-0	2	75	Undyed	36	Surgalloy™	UC-204
<b>CV-25</b> 1/2 Circle    22 mm 	4-0	1.5	75	Undyed	36	Surgalloy™	UC-403
	3-0	2	5 × 45	Undyed	12	D-Tach™	GC-33-M
	3-0	2	75	Undyed	36	Surgalloy™	UC-404
<b>GS-23</b> 1/2 Circle    22 mm 	2-0	3	75	Undyed	36		UC-863




## TAPER POINT

Needles Specifications	Size		Length cm	Color	Units Per box	Comments	Reorder code
	USP	EP					
<b>V-20</b> 1/2 Circle    26 mm 	4-0	1.5	75	Violet	36	Surgalloy™	GC-321
	4-0	1.5	75	Undyed	36	Surgalloy™	GC-121
	3-0	2	75	Undyed	36	Surgalloy™	GC-122
	3-0	2	75	Violet	36	Surgalloy™	GC-322
	2-0	3	75	Undyed	36	Surgalloy™	GC-123
	2-0	3	75	Violet	36	Surgalloy™	GC-323
	0	3.5	75	Undyed	36	Surgalloy™	GC-124
<b>GS-22</b> 1/2 Circle    27 mm 	2-0	3	75	Undyed	36		CC-883
	2-0	3	75	Violet	36		CC-830
	0	3.5	75	Undyed	36		CC-884
	0	3.5	75	Violet	36		CC-831
	1	4	75	Undyed	36		CC-885
<b>GS-21</b> 1/2 Circle    37 mm 	3-0	2	75	Undyed	36		CC-810
	3-0	2	90	Undyed	36		CC-922
	2-0	3	75	Violet	36		CC-843
	2-0	3	75	Undyed	36		CC-811
	2-0	3	90	Undyed	36		CC-923
	0	3.5	75	Violet	36		CC-844
	0	3.5	75	Undyed	36		CC-812
	0	3.5	90	Undyed	36		CC-924
	1	4	75	Violet	36		CC-845
	1	4	75	Undyed	36		CC-813
	1	4	90	Undyed	36		CC-925

## TAPER POINT

Needles Specifications	Size		Length cm	Color	Units Per box	Comments	Reorder code
	USP	EP					
<b>GS-24</b> 1/2 Circle    40 mm 	2-0	3	75	Undyed	36		CC-801
	2-0	3	90	Undyed	36		CC-915
	0	3.5	75	Undyed	36		CC-802
	0	3.5	90	Undyed	36		CC-913
	1	4	75	Undyed	36		CC-803
	1	4	90	Undyed	36		CC-914
<b>GS-25</b> 1/2 Circle    48 mm 	2-0	3	75	Undyed	36		CC-872
	0	3.5	90	Undyed	36		CC-904
	1	4	90	Violet	36		CC-975
	1	4	90	Undyed	36		CC-905
<b>GU-46</b> 5/8 Circle    27 mm 	3-0	2	75	Violet	36		UC-879
	2-0	3	75	Undyed	36		UC-978
<b>GU-46</b> 5/8 Circle    27 mm 	2-0	3	75	Violet	36	Surgalloy™	CC-658
<b>V-20</b> 1/2 Circle    26 mm Double-arm 							

## TAPERCUTTING

Needles Specifications	Size		Length cm	Color	Units Per box	Comments	Reorder code
	USP	EP					
<b>KV-34</b> 1/2 Circle    37 mm 	2-0	3	90	Undyed	12		CC-1945G
	1	4	90	Undyed	36	Surgalloy™	CC-196




<b>Description</b>	Biosyn™ monofilament synthetic absorbable sutures are prepared from a synthetic polyester composed of glycolide, dioxanone, and trimethylene carbonate. Biosyn™ synthetic absorbable sutures are colored violet to increase visibility and are also available undyed.
<b>Indications</b>	Biosyn™ synthetic absorbable sutures are indicated for use in general soft tissue approximation and/or ligation including use in ophthalmic surgery, but not for use in cardiovascular or neurological surgery.
<b>Actions</b>	Biosyn™ synthetic absorbable sutures elicit a minimal acute inflammatory reaction in tissues, which is followed by gradual encapsulation of the suture by fibrous connective tissue. Progressive loss of tensile strength and eventual absorption of Biosyn™ synthetic absorbable sutures occurs by means of hydrolysis where the suture is broken down to glycolic acid, dioxanoic acid, propane diol and carbon dioxide which are subsequently absorbed and metabolized by the body.
<b>Tensile Strength</b>	Studies in animals indicate tensile strength averages for Biosyn™ synthetic absorbable sutures are approximately 75% of U.S.P and E.P. minimum knot strength at two weeks and approximately 40% at three weeks post-implant.
<b>Sterilization</b>	Ethylene oxide.
<b>Absorbtion</b>	Absorption begins as a loss of tensile strength without appreciable loss of mass. Absorption of Biosyn™ synthetic absorbable sutures is essentially complete between 90 and 110 days.
<b>Packaging</b>	Biosyn™ synthetic absorbable sutures are available undyed (natural) or dyed violet in U.S.P. and E.P. sizes 1 (4 Metric) through 6-0 (0.7 Metric). The sutures are supplied sterile, in pre-cut lengths or affixed to various needle types using both permanent and removal needle attachment techniques. Please refer to the individual package IFUs for complete instructions, indications, contraindications, warnings and precautions.





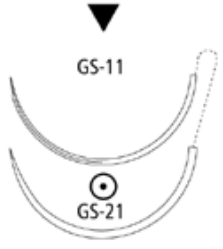

Biosyn™ suture is a brand you have trusted for over 15 years with a proven history of performance.

PRE-CUT							
Needles Specifications	Size		Length cm	Color	Units Per box	Comments	Reorder code
	USP	EP					
	<b>4-0</b>	1.5	6 × 45	Violet	24	Pre-cut	M-31-
	<b>3-0</b>	2	3 × 45	Violet	36	Pre-cut	M-72
	<b>3-0</b>	2	6 × 45	Violet	24	Pre-cut	M-32-
	<b>2-0</b>	3	6 × 45	Violet	24	Pre-cut	M-33-
	<b>2-0</b>	3	6 × 45	Undyed	24	Pre-cut	M-23-
	<b>0</b>	3.5	3 × 45	Violet	36	Pre-cut	M-74
	<b>1</b>	4	3 × 45	Violet	36	Pre-cut	M-75


STANDARD LENGHT							
Needles Specifications	Size		Length cm	Color	Units Per box	Comments	Reorder code
	USP	EP					
	<b>0</b>	3.5	150	Violet	36	Standard length	M-114

REVERSE CUTTING							
Needles Specifications	Size		Length cm	Color	Units Per box	Comments	Reorder code
	USP	EP					
<b>C-13</b> 3/8 Circle 19 mm 	<b>5-0</b>	1	75	Undyed	36		SM-690
	<b>4-0</b>	1.5	75	Violet	36		SM-791
	<b>4-0</b>	1.5	75	Undyed	12		SM-691G
	<b>4-0</b>	1.5	75	Undyed	36		SM-691
	<b>3-0</b>	2	75	Undyed	36		SM-822
	<b>3-0</b>	2	75	Violet	36		SM-922




## REVERSE CUTTING

Needles Specifications	Size USP	EP	Length cm	Color	Units Per box	Comments	Reorder code
<b>C-14</b> 3/8 Circle 24 mm 	<b>3-0</b>	2	75	Undyed	36		SM-693
	<b>2-0</b>	3	75	Violet	36		SM-923
	<b>2-0</b>	3	75	Undyed	36		SM-823
<b>C-23</b> 1/2 Circle 24 mm 	<b>4-0</b>	1.5	75	Violet	36		SM-633
	<b>3-0</b>	2	5 × 45	Violet	12	D-Tach™	SM-83-M
	<b>3-0</b>	2	75	Violet	36		SM-632
<b>GS-10</b> 1/2 Circle 26 mm 	<b>2-0</b>	3	75	Violet	36		CM-879
<b>GS-11</b> 1/2 Circle 37 mm 	<b>2-0</b>	3	75	Violet	36		CM-816
	<b>2-0</b>	3	75	Undyed	36		CM-543
	<b>0</b>	3.5	75	Violet	36		CM-817
	<b>0</b>	3.5	75	Undyed	36		CM-544
<b>GS-11 - GS-21</b> 1/2 Circle 37 mm *Double-arm 	<b>3-0</b>	2	135	Undyed	36		CM-991
	<b>2-0</b>	3	150	Undyed	36		CM-932
<b>HOS-11</b> 1/2 Circle 37 mm 	<b>1</b>	4	75	Violet	36		CM-528

## REVERSE CUTTING


Needles Specifications	Size		Length cm	Color	Units Per box	Comments	Reorder code
	USP	EP					
<b>GS-13</b> 1/2 Circle    48 mm 	<b>1</b>	4	90	Violet	36		CM-222

## PREMIUM REVERSE CUTTING


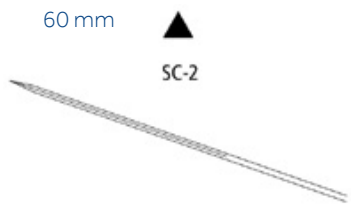
Needles Specifications	Size		Length cm	Color	Units Per box	Comments	Reorder code
	USP	EP					
<b>P-10</b> 3/8 Circle    11 mm 	<b>6-0</b>	0.7	45	Undyed	12	Surgalloy™	SM-3625
	<b>5-0</b>	1	45	Undyed	12	Surgalloy™	SM-3626
<b>P-13</b> 3/8 Circle    13 mm 	<b>6-0</b>	0.7	45	Undyed	12	Surgalloy™	SM-5688
	<b>5-0</b>	1	45	Undyed	12	Surgalloy™	SM-5687
	<b>4-0</b>	1.5	45	Undyed	12	Surgalloy™	SM-5690
<b>P-12</b> 3/8 Circle    19 mm 	<b>5-0</b>	1	45	Undyed	12	Surgalloy™	SM-5624
	<b>4-0</b>	1.5	45	Undyed	12	Surgalloy™	SM-5627
	<b>4-0</b>	1.5	75	Undyed	12	Surgalloy™	SM-5637
	<b>3-0</b>	2	45	Undyed	12	Surgalloy™	SM-5628
	<b>3-0</b>	2	75	Undyed	12	Surgalloy™	SM-5638






## PREMIUM REVERSE CUTTING

Needles Specifications	Size		Length cm	Color	Units Per box	Comments	Reorder code
	USP	EP					
<b>P-14</b> 3/8 Circle    24 mm 	4-0	1.5	45	Undyed	24	Surgalloy™	SM-5679
	4-0	1.5	75	Undyed	36	Surgalloy™	SM-5678
	3-0	2	75	Undyed	12	Surgalloy™	SM-5679G

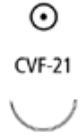

## CONVENTIONAL CUTTING

Needles Specifications	Size		Length cm	Color	Units Per box	Comments	Reorder code
	USP	EP					
<b>PC-13</b> 3/8 Circle    11 mm 	5-0	1	45	Undyed	24		CM-401
	4-0	1.5	45	Undyed	24		CM-400
<b>SC-2</b> Straight    60 mm 	4-0	1.5	75	Undyed	36		SM-644
	3-0	2	75	Undyed	36		SM-643

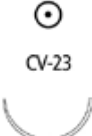


## X CUTTING

Needles Specifications	Size USP	EP	Length cm	Color	Units Per box	Comments	Reorder code
<b>DX-13</b> 3/8 Circle 13 mm 	5-0	1	45	Undyed	12	DermaX	XM-834
	4-0	1.5	45	Undyed	12	DermaX	XM-845
<b>DX-16</b> 3/8 Circle 16 mm 	5-0	1	45	Undyed	12	DermaX	XM-822
	4-0	1.5	45	Undyed	12	DermaX	XM-823
<b>DX-19</b> 3/8 Circle 19 mm 	5-0	1	45	Undyed	12	DermaX	XM-822
	4-0	1.5	45	Undyed	12	DermaX	XM-823





## TAPER POINT

Needles Specifications	Size USP	EP	Length cm	Color	Units Per box	Comments	Reorder code
<b>CVF-21</b> 1/2 Circle 12 mm 	5-0	1	45	Undyed	36	Surgalloy™	SM-432
	4-0	1.5	45	Undyed	36	Surgalloy™	SM-431
	4-0	1.5	75	Undyed	36	Surgalloy™	SM-434
<b>CV-11</b> 3/8 Circle 13 mm Double-arm 	5-0	1	75	Violet	36	Surgalloy™	CM-126



## TAPER POINT

Needles Specifications	Size USP	EP	Length cm	Color	Units Per box	Comments	Reorder code
<b>CV-22</b> 1/2 Circle 13 mm 	6-0	0.7	75	Violet	36	Surgalloy™	UM-986
	6-0	0.7	75	Undyed	36	Surgalloy™	UM-886
	5-0	1	75	Violet	36	Surgalloy™	GM-875
<b>CV-23</b> 1/2 Circle 17 mm 	6-0	0.7	75	Violet	36	Surgalloy™	UM-201
	5-0	1	75	Violet	36	Surgalloy™	UM-202
	5-0	1	75	Undyed	36	Surgalloy™	UM-213
	4-0	1.5	75	Undyed	36	Surgalloy™	UM-214
	4-0	1.5	75	Violet	36	Surgalloy™	UM-203
	3-0	2	75	Undyed	36	Surgalloy™	UM-215
	3-0	2	75	Violet	36	Surgalloy™	UM-204
	<b>CV-23</b> 1/2 Circle 17 mm Double-arm 	4-0	1.5	75	Violet	36	Surgalloy™
3-0		2	75	Violet	36	Surgalloy™	CM-208
<b>CV-25</b> 1/2 Circle 22 mm 	4-0	1.5	5 × 45	Violet	12	D-Tach™	GM-34-M
	4-0	1.5	75	Violet	36	Surgalloy™	GM-181
	3-0	2	45	Violet	12	D-Tach™	GMJ33M
	3-0	2	75	Violet	36	Surgalloy™	GM-324
	2-0	3	75	Violet	36	Surgalloy™	GM-183
<b>CV-25</b> 1/2 Circle 22 mm Double-arm 	4-0	1.5	75	Violet	36	Surgalloy™	GM-281
	4-0	1.5	90	Violet	36	Surgalloy™	GM-291
	3-0	2	90	Violet	36	Surgalloy™	GM-282






## TAPER POINT

Needles Specifications	Size USP	EP	Length cm	Color	Units Per box	Comments	Reorder code
<b>V-20</b> 1/2 Circle 26 mm 	4-0	1.5	75	Violet	36	Surgalloy™	GM-121
	4-0	1.5	75	Undyed	36	Surgalloy™	GM-321
	3-0	2	5 × 45	Undyed	12	D-Tach™	GM-69-M
	3-0	2	75	Violet	36	Surgalloy™	GM-122
	3-0	2	5 × 75	Violet	12	D-Tach™	GM-67-M
	3-0	2	75	Undyed	36	Surgalloy™	GM-322
	2-0	3	75	Violet	36		GM-123
	2-0	3	75	Undyed	36		GM-323
	2-0	3	90	Violet	36		GM-125
<b>V-20</b> 1/2 Circle 26 mm Double-arm 	4-0	1.5	75	Violet	36	Surgalloy™	GM-421
	4-0	1.5	90	Violet	36	Surgalloy™	GM-432
	3-0	2	90	Violet	36	Surgalloy™	GM-332
<b>GS-22</b> 1/2 Circle 27 mm 	3-0	2	75	Violet	36		CM-882
	2-0	3	75	Violet	36		CM-883
	2-0	3	75	Undyed	36		CM-830
	0	3.5	75	Violet	36		CM-884
	0	3.5	75	Undyed	36		CM-831
	1	4	75	Violet	36		CM-832
<b>V-30</b> 1/2 Circle 30 mm 	3-0	2	75	Violet	36		GM-222
	2-0	3	75	Violet	36		GM-223
	2-0	3	90	Violet	36		GM-326
	0	3.5	90	Violet	36		GM-327


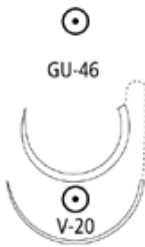

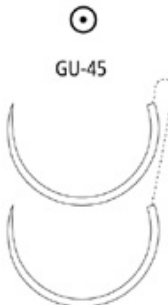

## TAPER POINT

Needles Specifications	Size		Length cm	Color	Units Per box	Comments	Reorder code
	USP	EP					
<b>V-26</b> 1/2 Circle 37 mm 	<b>3-0</b>	2	75	Undyed	36		GM-226
	<b>2-0</b>	3	75	Undyed	36		GM-227
<b>GS-21</b> 1/2 Circle 37 mm 	<b>3-0</b>	2	5 × 45	Violet	12	D-Tach™	CM-738-M
	<b>3-0</b>	2	75	Undyed	36		CM-842
	<b>3-0</b>	2	75	Violet	36		CM-810
	<b>3-0</b>	2	90	Violet	36		CM-922
	<b>2-0</b>	3	5 × 45	Violet	12	D-Tach™	CM-9M
	<b>2-0</b>	3	75	Undyed	36		CM-843
	<b>2-0</b>	3	75	Violet	36		CM-811
	<b>2-0</b>	3	90	Undyed	36		CM-953
	<b>2-0</b>	3	90	Violet	36		CM-923
	<b>0</b>	3.5	5 × 45	Undyed	12		CM-14-M
	<b>0</b>	3.5	5 × 45	Violet	12	D-Tach™	CM-10-M
	<b>0</b>	3.5	75	Undyed	36		CM-844
	<b>0</b>	3.5	75	Violet	36		CM-812
	<b>0</b>	3.5	90	Violet	36		CM-924
	<b>0</b>	3.5	90	Undyed	36		CM-954
	<b>1</b>	4	5 × 45	Violet	12	D-Tach™	CM-11-M
	<b>1</b>	4	75	Undyed	36		CM-845
	<b>1</b>	4	75	Violet	36		CM-813
	<b>1</b>	4	90	Violet	36		CM-925



## TAPER POINT

Needles Specifications	Size		Length cm	Color	Units Per box	Comments	Reorder code
	USP	EP					
<b>GS-24</b> 1/2 Circle 40 mm 	<b>3-0</b>	2	90	Undyed	36		CM-962
	<b>2-0</b>	3	75	Violet	36		CM-801
	<b>0</b>	3.5	90	Violet	36		CM-914
	<b>0</b>	3.5	90	Undyed	36		CM-964
	<b>1</b>	4	90	Violet	36		CM-966
<b>GS-24 LOOP</b> 1/2 Circle 40 mm 	<b>0</b>	3.5	150	Violet	12	Loop	CM-765-L
	<b>1</b>	4	150	Violet	12	Loop	CM-767-L
<b>GS-25</b> 1/2 Circle 48 mm 	<b>2-0</b>	3	90	Violet	36		CM-903
	<b>0</b>	3.5	90	Undyed	36		CM-974
	<b>0</b>	3.5	90	Violet	36		CM-904
	<b>1</b>	4	5 × 45	Violet	12	D-Tach™	CM-31-M
	<b>1</b>	4	90	Undyed	36		CM-975
	<b>1</b>	4	90	Violet	36		CM-905
<b>GS-25 LOOP</b> 1/2 Circle 48 mm 	<b>0</b>	3.5	150	Violet	12	Loop	CM-790L
	<b>1</b>	4	150	Violet	12	Loop	CM-776L
<b>GS-26 LOOP</b> 1/2 Circle 65 mm 	<b>0</b>	3.5	150	Violet	12	Loop	CM-779L
	<b>1</b>	4	240	Violet	12	Loop	CM-780L

## TAPER POINT



Needles Specifications	Size USP	EP	Length cm	Color	Units Per box	Comments	Reorder code
<b>GU-46</b> 5/8 Circle 27 mm 	3-0	2	75	Violet	36		UM-879
	2-0	3	75	Undyed	36		UM-978
	2-0	3	75	Violet	36		UM-878
	0	3.5	75	Violet	36		UM-877
	0	3.5	75	Undyed	36		UM-977
<b>GU-46 - V - 20</b> 5/8 Circle 27 mm 	2-0	3	75	Violet	36	Loop	CM-558
<b>GU-45</b> 5/8 Circle 37 mm 	2-0	3	75	Undyed	36		UM-345
	2-0	3	75	Violet	36		UM-245
<b>GU-45 - GU-45</b> 5/8 Circle 37 mm Double-arm (Not to scale) 	3-0	2	75	Violet	36		UM-243
<b>GU-44</b> 5/8 Circle 40 mm 	0	3.5	75	Violet	36		UM-383

## HEAVY TAPER POINT

Needles Specifications	Size		Length cm	Color	Units Per box	Comments	Reorder code
	USP	EP					
<b>HGS-20</b> 1/2 Circle    34 mm 	<b>0</b>	3.5	5 × 45	Violet	12	D-Tach™	CM-90-M
<b>HGS-21</b> 1/2 Circle    37 mm 	<b>0</b>	3.5	75	Violet	36		CM-436



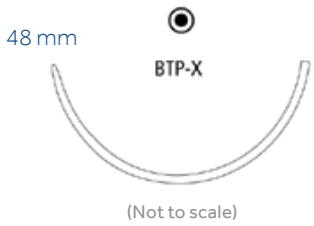
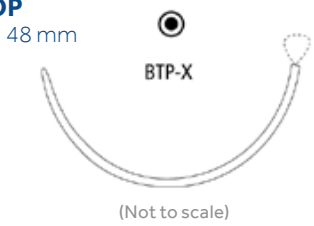
## TAPERCUTTING

Needles Specifications	Size		Length cm	Color	Units Per box	Comments	Reorder code
	USP	EP					
<b>KV-34</b> 1/2 Circle 37 mm 	<b>2-0</b>	3	75	Undyed	36		CM-827
<b>KV-40 LOOP</b> 1/2 Circle 48 mm 	<b>1</b>	4	150	Violet	12	Loop	CM-740L

## BLUNT POINT

Needles Specifications	Size		Length cm	Color	Units Per box	Comments	Reorder code
	USP	EP					
<b>BTV-20</b> 1/2 Circle 26 mm 	<b>3-0</b>	2	75	Undyed	36		GMT-226
<b>BTP-1</b> 1/2 Circle 37 mm 	<b>2-0</b>	3	75	Violet	36	ProtectPoint™	CMT-811
	<b>0</b>	3.5	5 × 45	Violet	12	D-Tach™, ProtectPoint™	CMT-740-M
	<b>0</b>	3.5	75	Violet	36		CMT-812
	<b>0</b>	3.5	90	Undyed	36	ProtectPoint™	CMT-967
	<b>1</b>	4	5 × 45	Violet	12	D-Tach™, ProtectPoint™	CMT-741-M
<b>BGS-24</b> 1/2 Circle 40 mm 	<b>1</b>	4	90	Violet	36	ProtectPoint™	CM-951
<b>BGS-25</b> 1/2 Circle 48 mm 	<b>1</b>	4	5 × 45	Violet	12	D-Tach™, ProtectPoint™	CM-949-M

## BLUNT POINT

Needles Specifications	Size		Length cm	Color	Units Per box	Comments	Reorder code
	USP	EP					
<b>BTP-X</b> 1/2 Circle 48 mm  (Not to scale)	<b>0</b>	3.5	5 × 45	Violet	12	D-Tach™, ProtectPoint™	CMT-950-M
	<b>0</b>	3.5	90	Undyed	36	ProtectPoint™	CMT-974
	<b>1</b>	4	5 × 45	Violet	12	D-Tach™, ProtectPoint™	CMT-949-M
<b>BTP-X LOOP</b> 1/2 Circle 48 mm  (Not to scale)	<b>0</b>	3.5	150	Violet	12	LOOP, ProtectPoint™	CMT-789-L

# MAXON™

## Monofilament Absorbable Suture

Highest knot holding capacity versus competition.






Please refer to the individual package IFUs for complete instructions, indications, contraindications, warnings and precautions.






<b>Description</b>	Maxon™ monofilament polyglyconate synthetic absorbable sutures are prepared from a copolymer of glycolic acid and trimethylene carbonate.
<b>Indications</b>	Maxon™ sutures are indicated for use as absorbable sutures in general soft tissue approximation and/or ligation, including use in pediatric cardiovascular tissue, where growth is expected to occur, and in peripheral vascular tissue. Maxon™ sutures are not indicated for use in adult cardiovascular tissue, ophthalmic surgery, microsurgery, and neural tissue.
<b>Actions</b>	Maxon™ synthetic absorbable sutures elicit a minimal acute inflammatory reaction in tissues, which is followed by gradual encapsulation of the suture by fibrous connective tissue. Progressive loss of tensile strength and eventual absorption of Maxon™ synthetic absorbable sutures occurs by means of hydrolysis where the suture is broken down and subsequently absorbed and metabolized by the body.
<b>Tensile Strength</b>	Studies indicate that tensile strength for Maxon™ sutures is approximately 80% of the original suture strength at one week post-implant. Approximately 75% of the original suture strength remains at two weeks, 65% at three weeks, 50% at four weeks and 25% at six weeks post-implant.
<b>Sterilization</b>	Ethylene oxide.
<b>Absorbtion</b>	Absorption of Maxon™ synthetic absorbable sutures is minimal until about the 60th day after implantation and is essentially complete within six months.
<b>Packaging</b>	Maxon™ sutures are available in USP and E.P sizes 1 (4 Metric) through 7-0 (0.5 Metric). Maxon™ sutures are dyed with D&C Green No. 6. Maxon™ sutures are also available in clear (Undyed). The sutures are supplied sterile, affixed to various needle types, using permanent needle attachment techniques. Please refer to the individual package IFUs for complete instructions, indications, contraindications, warnings and precautions.

Maxon™ suture is a brand you have trusted for over 25 years with a history of improvements and performance





PRE-CUT							
Needles Specifications	Size		Length cm	Color	Units Per box	Comments	Reorder code
	USP	EP					
	<b>3-0</b>	2	3 × 45	Green	36	Pre-cut	8886600741
	<b>2-0</b>	3	3 × 45	Green	36	Pre-cut	8886600751
	<b>0</b>	3.5	5 × 45	Green	12	Pre-cut	8886601563

REVERSE CUTTING							
Needles Specifications	Size		Length cm	Color	Units Per box	Comments	Reorder code
	USP	EP					
<b>C-1</b> 3/8 Circle 12 mm 	<b>5-0</b>	1	45	Green	36		8886660521
	<b>4-0</b>	1.5	45	Green	36		8886660531
<b>C-12</b> 3/8 Circle 16 mm 	<b>4-0</b>	1.5	75	Clear	36		8886611831
<b>C-13</b> 3/8 Circle 19 mm 	<b>5-0</b>	1	45	Clear	36		8886660821
	<b>4-0</b>	1.5	45	Clear	36		8886660831
	<b>4-0</b>	1.5	75	Green	36		8886663931
	<b>4-0</b>	1.5	75	Clear	36		8886661231
	<b>3-0</b>	2	45	Clear	36		8886660841
	<b>3-0</b>	2	75	Green	36		8886663941
	<b>3-0</b>	2	75	Clear	36		8886661241

## REVERSE CUTTING






Needles Specifications	Size		Length cm	Color	Units Per box	Comments	Reorder code
	USP	EP					
<b>C-14</b> 3/8 Circle 24 mm 	4-0	1.5	75	Clear	36		8886661831
	3-0	2	75	Clear	36		8886661841
	3-0	2	75	Green	36		8886664141
	2-0	3	75	Green	36		8886664151
	2-0	3	75	Clear	36		8886661851
<b>C-16</b> 3/8 Circle 30 mm 	2-0	3	75	Clear	36		8886662651
<b>C-27</b> 1/2 Circle 20 mm 	3-0	2	45	Green	36		8886605941
<b>GS-10</b> 1/2 Circle 26 mm 	3-0	2	75	Clear	36		8886674441
	2-0	3	5 × 45	Green	12	D-Tach™	CMM-435-MG
	2-0	3	75	Green	36		8886674551
	2-0	3	75	Clear	36		8886674451
<b>HOS-10</b> 1/2 Circle 26 mm 	1	4	75	Green	36		88866751-71

## REVERSE CUTTING

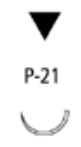

Needles Specifications	Size		Length cm	Color	Units Per box	Comments	Reorder code
	USP	EP					
<b>GS-11</b> 1/2 Circle 37 mm 	<b>2-0</b>	3	75	Green	36		8886675551
	<b>0</b>	3.5	75	Green	36		8886675561
	<b>1</b>	4	75	Green	36		8886675571
<b>HOS-11</b> 1/2 Circle 37 mm 	<b>1</b>	4	75	Green	36		88866757-71
<b>HOS-12</b> 1/2 Circle 40 mm 	<b>1</b>	4	75	Green	36		8886678771
<b>GS-13</b> 1/2 Circle 48 mm 	<b>1</b>	4	90	Green	36		8886676971



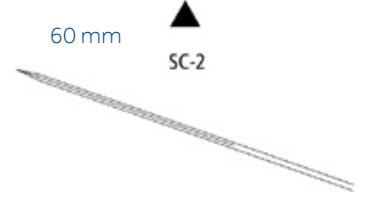
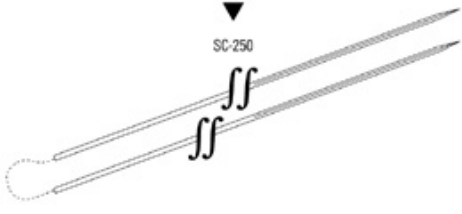
## PREMIUM REVERSE CUTTING

Needles Specifications	Size		Length cm	Color	Units Per box	Comments	Reorder code
	USP	EP					
<b>P-10</b> 3/8 Circle    11 mm 	6-0	0.7	45	Green	36	Surgalloy™	8886653311
	6-0	0.7	45	Clear	36	Surgalloy™	8886653211
	5-0	1	45	Clear	36	Surgalloy™	8886653221
<b>P-13</b> 3/8 Circle    13 mm 	6-0	0.7	45	Clear	36	Surgalloy™	SMM-5516
	5-0	1	45	Clear	36	Surgalloy™	SMM-5526
	4-0	1.5	45	Clear	36	Surgalloy™	SMM-5536
<b>P-11</b> 3/8 Circle    16 mm 	5-0	1	45	Clear	36	Surgalloy™	SMM-5221
	4-0	1.5	45	Clear	36	Surgalloy™	SMM-5431
<b>P-12</b> 3/8 Circle    19 mm 	5-0	1	45	Clear	36	Surgalloy™	SMM-5042
	4-0	1.5	45	Clear	36	Surgalloy™	SMM-5043
	4-0	1.5	45	Green	36	Surgalloy™	SMM-5033
	4-0	1.5	75	Clear	36	Surgalloy™	SMM-5323
	3-0	2	45	Clear	36	Surgalloy™	SMM-5044
	3-0	2	75	Clear	36	Surgalloy™	8886663241
<b>P-14</b> 3/8 Circle    24 mm 	3-0	2	45	Clear	36	Surgalloy™	SMM-5241
	2-0	3	45	Clear	36		8886662451


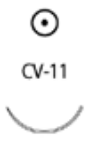
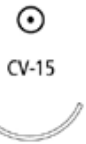

## PREMIUM REVERSE CUTTING

Needles Specifications	Size		Length cm	Color	Units Per box	Comments	Reorder code
	USP	EP					
<b>P-21</b> 1/2 Circle 9 mm 	5-0	1	45	Clear	12		88866704-23
<b>P-24</b> 1/2 Circle 16 mm 	4-0	1.5	45	Clear	36	Surgalloy™	SMM-5831

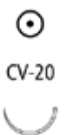




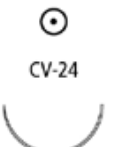
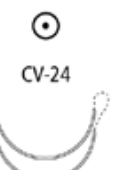
## CONVENTIONAL CUTTING

Needles Specifications	Size		Length cm	Color	Units Per box	Comments	Reorder code
	USP	EP					
<b>SC-2</b> Straight 60 mm 	4-0	1.5	75	Clear	36		8886690831
	3-0	2	75	Clear	36		8886690841
<b>SC-250</b> Straight 254 mm Double-arm (Not to scale) 	2-0	3	75	Green	12		8886639053



## TAPER POINT

Needles Specifications	Size USP	EP	Length cm	Color	Units Per box	Comments	Reorder code	
<b>CV-1</b> 3/8 Circle 9 mm Double-arm		7-0	0.5	45	Green	36	Surgalloy™	8886614501
		7-0	0.5	60	Green	36	Surgalloy™	8886613901
		7-0	0.5	75	Green	36	Surgalloy™	8886612701
		6-0	0.7	60	Green	36	Surgalloy™	8886613911
		6-0	0.7	75	Green	36	Surgalloy™	8886612711
<b>CV</b> 3/8 Circle 11 mm Double-arm		6-0	0.7	45	Green	36	Surgalloy™	8886611511
<b>CV-11</b> 3/8 Circle 13 mm		6-0	0.7	75	Green	36	Surgalloy™	8886614811
		5-0	1	75	Green	36	Surgalloy™	8886614821
<b>CV-11</b> 3/8 Circle 13 mm Double-arm		6-0	0.7	45	Green	36	Surgalloy™	8886612111
		6-0	0.7	60	Green	36	Surgalloy™	8886614711
		6-0	0.7	75	Green	36	Surgalloy™	8886612411
		5-0	1	4×60	Green	12	Surgalloy™	8886616523
		5-0	1	60	Green	36	Surgalloy™	8886614721
		5-0	1	75	Green	36	Surgalloy™	8886612421
		5-0	1	90	Green	36	Surgalloy™	8886612821
<b>CV-15</b> 3/8 Circle 17 mm		4-0	1.5	75	Green	36		8886605731
<b>CV-15</b> 3/8 Circle 17 mm Double-arm		5-0	1	75	Green	36	Surgalloy™	8886646221
		4-0	1.5	75	Green	36	Surgalloy™	8886646231

## TAPER POINT

Needles Specifications	Size USP	EP	Length cm	Color	Units Per box	Comments	Reorder code
<b>CV-20</b> 1/2 Circle 10 mm 	6-0	0.7	75	Green	36	Surgalloy™	8886600111
<b>CV-22</b> 1/2 Circle 13 mm 	6-0	0.7	75	Green	36	Surgalloy™	8886621311
	5-0	1	75	Green	36		8886621321
	4-0	1.5	5 × 45	Green	12	D-Tach™	8886620133
	4-0	1.5	75	Green	36	Surgalloy™	8886621331
<b>CV-22</b> 1/2 Circle 13 mm Double-arm 	6-0	0.7	75	Green	36	Surgalloy™	8886613211
	5-0	1	75	Green	36	Surgalloy™	8886612221
<b>CV-23</b> 1/2 Circle 17 mm 	5-0	1	75	Green	36	Surgalloy™	8886620721
	4-0	1.5	75	Green	36	Surgalloy™	8886620731
	3-0	2	75	Green	36	Surgalloy™	8886620741
<b>CV-23</b> 1/2 Circle 17 mm Double-arm 	5-0	1	75	Green	36	Surgalloy™	8886614921
	4-0	1.5	90	Green	36	Surgalloy™	8886613431
<b>CV-24</b> 1/2 Circle 20 mm 	5-0	1	75	Green	36		8886621721
	4-0	1.5	5 × 45	Green	12	D-Tach™	8886621933
	4-0	1.5	75	Green	36		8886621731
	3-0	2	5 × 45	Green	12	D-Tach™	8886621943
	3-0	2	75	Green	36		8886621741
<b>CV-24</b> 1/2 Circle 20 mm Double-arm 	4-0	1.5	90	Green	36		8886622131
	3-0	2	90	Green	36		8886603641



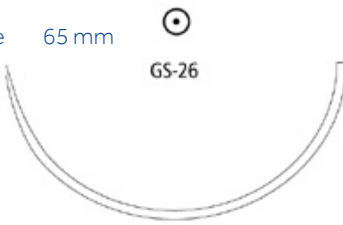



## TAPER POINT

Needles Specifications	Size		Length cm	Color	Units Per box	Comments	Reorder code
	USP	EP					
<b>V-20</b> 1/2 Circle 26 mm 	4-0	1.5	5 × 45	Green	12	D-Tach™	8886622933
	4-0	1.5	5 × 75	Green	12	D-Tach™	8886603433
	4-0	1.5	75	Green	36	Surgalloy™	8886623331
	3-0	2	5 × 45	Green	12	D-Tach™	8886622943
	3-0	2	75	Green	36	Surgalloy™	8886623341
	3-0	2	75	Clear	36	Surgalloy™	8886623241
	3-0	2	5 × 75	Green	12	D-Tach™	8886603443
	2-0	3	5 × 45	Green	12	D-Tach™	8886622953
	2-0	3	75	Green	36		8886623351
<b>V-20</b> 1/2 Circle 26 mm Double-arm 	4-0	1.5	75	Green	36	Surgalloy™	8886623131
	4-0	1.5	90	Green	36		8886623531
	4-0	1.5	105	Green	36		8886614331
	3-0	2	90	Green	36		8886623541
<b>GS-22</b> 1/2 Circle 27 mm 	3-0	2	75	Green	36		8886624541
	2-0	3	5 × 45	Green	12	D-Tach™	8886624753
	2-0	3	75	Green	36		8886624551
	0	3.5	75	Green	36		8886624561
	1	4	75	Green	36		8886624571
<b>V-30</b> 1/2 Circle 30 mm 	4-0	1.5	75	Green	36		8886626131
	3-0	2	75	Green	36		8886626141
	2-0	3	75	Green	36		8886626151
	0	3.5	75	Green	36		8886626161

## TAPER POINT

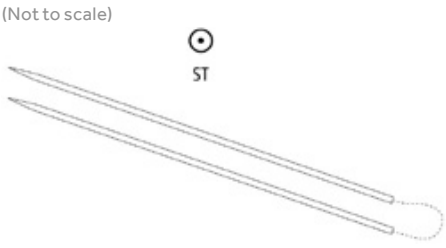
Needles Specifications	Size		Length cm	Color	Units Per box	Comments	Reorder code	
	USP	EP						
<b>GS-21</b> 1/2 Circle 37 mm 	3-0	2	75	Green	36		8886626741	
	2-0	3	5 × 45	Green	12	D-Tach™	8886627553	
	2-0	3	75	Green	36		8886626751	
	2-0	3	90	Green	36		8886626951	
	2-0	3	150	Green	12		GMM-951	
	0	3.5	5 × 45	Green	12	D-Tach™	8886627563	
	0	3.5	75	Green	36		8886626761	
	0	3.5	90	Green	36		8886626961	
	1	4	5 × 45	Green	12	D-Tach™	8886627573	
	1	4	75	Green	36		8886626771	
	1	4	90	Green	36		8886626971	
	<b>V-26</b> 1/2 Circle 37 mm 	3-0	2	75	Green	36		8886625741
2-0		3	75	Green	36		8886625751	
0		3.5	75	Green	36		8886625761	
<b>GS-24</b> 1/2 Circle 40 mm 	2-0	3	90	Green	36		8886628751	
	0	3.5	75	Green	36		8886628561	
	0	3.5	90	Green	36		8886628761	
	1	4	75	Green	36		8886628571	
	1	4	90	Green	36		8886628771	
<b>GS-24 LOOP</b> 1/2 Circle 40 mm 	0	3.5	150	Green	12	Loop	GMM-563L	
	1	4	150	Green	12	Loop	GMM573L	

## TAPER POINT

Needles Specifications	Size		Length cm	Color	Units Per box	Comments	Reorder code
	USP	EP					
<b>GS-25</b> 1/2 Circle 48 mm 	0	3.5	5 × 45	Green	12	D-Tach™	8886631163
	0	3.5	90	Green	36		8886630761
	1	4	5 × 45	Green	12	D-Tach™	8886631173
	1	4	90	Green	36		8886630771
	1	4	150	Green	24		8886630972
<b>GS-25 LOOP</b> 1/2 Circle 48 mm 	0	3.5	150	Green	12	Loop	GMM-634L
	1	4	150	Green	12	Loop	GMM-632L
<b>GS-26</b> 1/2 Circle 65 mm 	0	3.5	150	Green	24		8886631962
	1	4	100	Green	36		8886632271
	1	4	150	Green	24		8886631972
<b>GS-26 LOOP</b> 1/2 Circle 65 mm 	0	3.5	150	Green	12	Loop	GMM-341L
	0	3.5	240	Green	12	Loop	GMM-343L
	1	4	240	Green	12	Loop	GMM-344L
<b>GU-46</b> 5/8 Circle 27 mm 	2-0	3	75	Green	36		8886640351
<b>GU-46</b> 5/8 Circle 27 mm Double-arm 	2-0	3	75	Green	36		8886640451



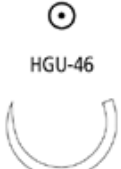
## TAPER POINT

Needles Specifications	Size		Length cm	Color	Units Per box	Comments	Reorder code
	USP	EP					
<b>ST</b> Straight 60 mm Double-arm (Not to scale)	3-0	2	75	Green	36		8886623741
	0	3.5	75	Green	36		8886623761







## HEAVY TAPER POINT

Needles Specifications	Size		Length cm	Color	Units Per box	Comments	Reorder code
	USP	EP					
<b>HGS-22</b> 1/2 Circle 27 mm	0	3.5	75	Green	36		8886624961
<b>HGS-21</b> 1/2 Circle 37 mm	0	3.5	90	Green	36		8886633961
	1	4	90	Green	36		88866339-71
<b>HGU-46</b> 5/8 Circle 27 mm	0	3.5	45	Green	36		8886645561
	0	3.5	75	Green	36		8886645661
	1	4	75	Green	36		8886641971

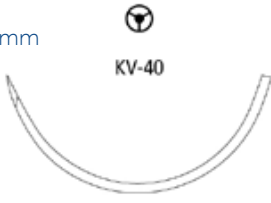








## TAPERCUTTING

Needles Specifications	Size USP	EP	Length cm	Color	Units Per box	Comments	Reorder code
<b>KV-1</b> 3/8 Circle    9 mm 	6-0	0.7	75	Green	36	Surgalloy™	88866768-11
<b>KV-11</b> 3/8 Circle    13 mm 	6-0	0.7	75	Green	36	Surgalloy™	88866297-11
	5-0	1	75	Green	36	Surgalloy™	88866297-21
<b>KV-11</b> 3/8 Circle    13 mm Double-arm 	6-0	0.7	75	Green	36	Surgalloy™	88866288-11
	5-0	1	75	Green	36	Surgalloy™	88866288-21
<b>KV-5</b> 1/2 Circle    17 mm Double-arm 	4-0	1.5	90	Green	36	Surgalloy™	VMM-137
<b>KV-7</b> 1/2 Circle    26 mm 	2-0	3	75	Green	36		88866292-51
<b>KV-34</b> 1/2 Circle    37 mm 	0	3.5	90	Green	36		8886627761





## TAPERCUTTING

Needles Specifications	Size		Length cm	Color	Units Per box	Comments	Reorder code
	USP	EP					
<b>KV-37 LOOP</b> 1/2 Circle 40 mm 	1	4	150	Green	12	Loop	GMM-873L
<b>KV-40</b> 1/2 Circle 48 mm 	1	4	90	Green	36		8886629571
<b>KV-40 LOOP</b> 1/2 Circle 48 mm 	0	3.5	150	Green	12	Loop	GMM-806L
	1	4	150	Green	12	Loop	GMM-807L

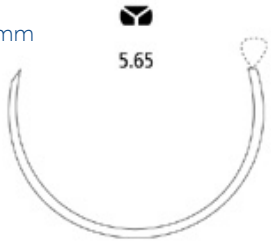
## BLUNT POINT

Needles Specifications	Size		Length cm	Color	Units Per box	Comments	Reorder code
	USP	EP					
<b>BTP-1</b> 1/2 Circle 37 mm 	2-0	3	90	Green	36	ProtectPoint™	GMMT-951
	0	3.5	5 × 45	Green	12	D-Tach™, ProtectPoint™	GMMT-563- MG
	0	3.5	90	Green	36	ProtectPoint™	GMMT-961


## BLUNT POINT

Needles Specifications	Size		Length cm	Color	Units Per box	Comments	Reorder code
	USP	EP					
<b>BTP-X</b> 1/2 Circle 48 mm 	0	3.5	5 × 45	Green	12	D-Tach™, ProtectPoint™	GMMT-541-MG
	1	4	5 × 45	Green	12	D-Tach™, ProtectPoint™	GMMT-540-MG
<b>BTP-X LOOP</b> 1/2 Circle 48 mm 	0	3.5	150	Green	12	Loop, ProtectPoint™	GMMT-89L
	1	4	150	Green	12	Loop, ProtectPoint™	GMMT-90L
<b>BGS-25</b> 1/2 Circle 48 mm 	0	3.5	90	Green	12	ProtectPoint™	8886606163
	1	4	90	Green	12	ProtectPoint™	8886606173
<b>BGS-29</b> 1/2 Circle 63 mm 	1	4	150	Green	12	ProtectPoint™	8886606773

## MULTIPURPOSE

Needles Specifications	Size		Length cm	Color	Units Per box	Comments	Reorder code
	USP	EP					
<b>5.65 LOOP</b> 5/8 Circle 65 mm 	0	3.5	150	Green	12	Loop	GMM-496L
	1	4	180	Green	12	Loop	GMM-497L

## SPATULA

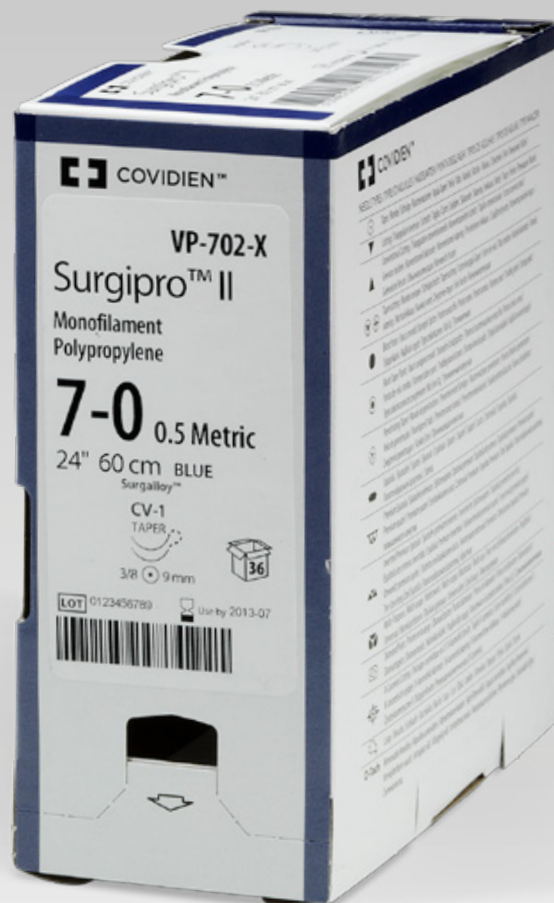
Needles Specifications	Size		Length cm	Color	Units Per box	Comments	Reorder code
	USP	EP					
<b>DGE-6</b> 3/8 Circle 8 mm 	7-0	0.5	45	Green	36		8886662801



# SURGIPRO™ II

Monofilament Polypropylene

Meeting the demands in cardiovascular surgery.



Please refer to the individual package IFUs for complete instructions, indications, contraindications, warnings and precautions.

**Description**

Surgipro™ and Surgipro™ II monofilament sutures (pigmented) are inert, nonabsorbable, sterile sutures composed of an isotactic crystalline stereoisomer of polypropylene (a synthetic linear polyolefin) and polyethylene.

**Indications**

Surgipro™ and Surgipro™ II polypropylene sutures are indicated for use in general soft tissue approximation and/or ligation, including use in cardiovascular, ophthalmic, and neural tissue.

**Actions**

Surgipro™ and Surgipro™ II polypropylene sutures elicit a minimal acute inflammatory reaction in tissues, which is followed by gradual encapsulation of the suture by fibrous connective tissue.

**Tensile Strength**

Surgipro™ and Surgipro™ II polypropylene sutures are not absorbed, nor is any significant change in strength retention known to occur in vivo.

**Sterilization**




Ethylene oxide.

**Packaging**

Surgipro™ II polypropylene sutures are available either undyed (clear) or dyed with Copper Phthalocyanine Blue. Surgipro™ sutures are available in sizes 1 (4 Metric), through 6-0 (0.7 Metric). Surgipro™ II sutures are available in sizes 3-0 (2 Metric) through 8-0 (0.4 Metric). The sutures are supplied sterile; affixed to various types of needles, using both permanent needle attachment techniques. Surgipro™ and Surgipro™ II polypropylene sutures are also available with PTFE (polytetrafluorethylene) Pledgets for use as a pad between the suture area and the tissue surface to increase the load bearing area. Surgipro™ polypropylene sutures are available in presentations containing Bead and Collar components to anchor the ends of the suture for subcuticular closure or for use as tendon sutures.




## Surgipro™ / Surgipro™ II

Surgipro™ / Surgipro™ II suture is a brand you have trusted for many years with several product innovations.

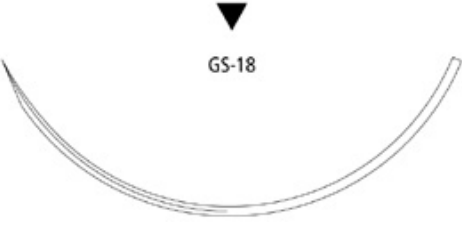



REVERSE CUTTING							
Needles Specifications	Size USP	EP	Length cm	Color	Units Per box	Comments	Reorder code
<b>C-12</b> 3/8 Circle 16 mm 	<b>6-0</b>	0.7	45	Blue	12		SP-660G
	<b>6-0</b>	0.7	45	Blue	36		SP-660
	<b>5-0</b>	1	45	Blue	12		SP-559G
	<b>4-0</b>	1.5	75	Blue	36		SP-631
<b>C-13</b> 3/8 Circle 19 mm 	<b>5-0</b>	1	45	Blue	36		SP-661
	<b>5-0</b>	1	45	Blue	12		SP-661G
	<b>4-0</b>	1.5	45	Blue	36		SP-683
	<b>4-0</b>	1.5	45	Blue	12		SP-683G
	<b>3-0</b>	2	45	Blue	12		SP-665G
	<b>3-0</b>	2	45	Blue	36		SP-665
	<b>3-0</b>	2	75	Blue	36		SP-655
<b>C-14</b> 3/8 Circle 24 mm 	<b>4-0</b>	1.5	45	Blue	36		SP-629
	<b>3-0</b>	2	45	Blue	36		SP-684
	<b>3-0</b>	2	45	Blue	12		SP-684G








## REVERSE CUTTING

Needles Specifications	Size		Length cm	Color	Units Per box	Comments	Reorder code
	USP	EP					
<b>C-15</b> 3/8 Circle    26 mm 	<b>2-0</b>	3	45	Blue	12		SP-868G
	<b>2-0</b>	3	45	Blue	36		SP-868
	<b>2-0</b>	3	75	Blue	36		SP-685
<b>C-16</b> 3/8 Circle    30 mm 	<b>4-0</b>	1.5	75	Blue	36	Nupack	SP-670
	<b>3-0</b>	2	75	Blue	36		SP-675
<b>C-18</b> 3/8 Circle    35 mm 	<b>3-0</b>	2	75	Blue	12	2 beads & 2 collars	SP-624-BC
<b>C-17</b> 3/8 Circle    39 mm 	<b>3-0</b>	2	75	Blue	36		SP-649
	<b>2-0</b>	3	75	Blue	36		SP-689
	<b>0</b>	3.5	75	Blue	36		SP-690




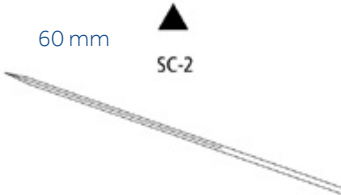

## REVERSE CUTTING

Needles Specifications	Size		Length cm	Color	Units Per box	Comments	Reorder code
	USP	EP					
<b>GS-18</b> 3/8 Circle 77 mm 	<b>2</b>	5	75	Blue	24		CP-460
<b>GS-18</b> 3/8 Circle 77 mm Double-arm 	<b>2</b>	5	3 × 60	Blue	12	Long pack	3-CP-360
<b>GS-10</b> 1/2 Circle 26 mm 	<b>2-0</b>	3	75	Blue	36		CP-443
<b>GS-11</b> 1/2 Circle 37 mm 	<b>2-0</b>	3	75	Blue	36		CP-453


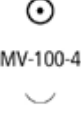


## PREMIUM REVERSE CUTTING

Needles Specifications	Size USP	EP	Length cm	Color	Units Per box	Comments	Reorder code
<b>P-10</b> 3/8 Circle 11 mm 	<b>7-0</b>	0.5	45	Blue	12	Surgalloy™	SP-1696G
	<b>7-0</b>	0.5	45	Blue	36	Surgalloy™	SP-1696
	<b>6-0</b>	0.7	45	Blue	36	Surgalloy™	SP-1697
	<b>6-0</b>	0.7	45	Clear	36	Surgalloy™	SP-1606
	<b>6-0</b>	0.7	45	Blue	12	Surgalloy™	SP-5697G
	<b>5-0</b>	1	45	Blue	12	Surgalloy™	SP-2698G
<b>P-13</b> 3/8 Circle 13 mm 	<b>6-0</b>	0.7	45	Blue	12	Surgalloy™	SP-5695-G
	<b>6-0</b>	0.7	45	Blue	36	Surgalloy™	SP-5695
	<b>5-0</b>	1	45	Blue	12	Surgalloy™	SP-5698-G
	<b>5-0</b>	1	45	Blue	36	Surgalloy™	SP-5698
	<b>4-0</b>	1.5	45	Blue	12	Surgalloy™	SP-5699G
	<b>4-0</b>	1.5	45	Blue	36	Surgalloy™	SP-5699
<b>P-11</b> 3/8 Circle 16 mm 	<b>6-0</b>	0.7	45	Blue	12	Surgalloy™	SP-5680G
	<b>5-0</b>	1	45	Blue	12	Surgalloy™	SP-5681G
	<b>5-0</b>	1	45	Blue	36	Surgalloy™	SP-5681
	<b>4-0</b>	1.5	45	Blue	12	Surgalloy™	SP-5692G
	<b>3-0</b>	2	45	Blue	12	Surgalloy™	SP-1693G
<b>P-12</b> 3/8 Circle 19 mm 	<b>5-0</b>	1	45	Blue	12	Surgalloy™	SP-5686G
	<b>5-0</b>	1	45	Blue	36	Surgalloy™	SP-5686
	<b>4-0</b>	1.5	45	Blue	12	Surgalloy™	SP-5682G
	<b>4-0</b>	1.5	45	Blue	36	Surgalloy™	SP-5682
	<b>3-0</b>	2	45	Blue	36	Surgalloy™	SP-5687
	<b>3-0</b>	2	45	Blue	12	Surgalloy™	SP-5687G
	<b>3-0</b>	2	75	Blue	36		SP-5587
<b>P-14</b> 3/8 Circle 24 mm 	<b>3-0</b>	2	45	Blue	12	Surgalloy™	SP-5663G
	<b>3-0</b>	2	75	Blue	36		SP-5563

## CONVENTIONAL CUTTING

Needles Specifications	Size		Length cm	Color	Units Per box	Comments	Reorder code
	USP	EP					
<b>PC-13</b> 3/8 Circle 11 mm 	<b>6-0</b>	0.7	45	Blue	36		SP-1617
	<b>5-0</b>	1	45	Blue	36		SP-1618
<b>PC-11</b> 3/8 Circle 16 mm 	<b>6-0</b>	0.7	45	Blue	36		SP-1636
	<b>5-0</b>	1	45	Blue	36	Surgalloy™	SP-1635
	<b>4-0</b>	1.5	45	Blue	36		SP-1634
<b>PC-12</b> 3/8 Circle 19 mm 	<b>4-0</b>	1.5	45	Blue	36		SP-1644
	<b>3-0</b>	2	45	Blue	36		SP-1643
<b>SC-2</b> Straight 60 mm 	<b>4-0</b>	1.5	75	Blue	36		SP-621
	<b>3-0</b>	2	75	Blue	36		SP-622
	<b>3-0</b>	2	75	Blue	24	2 beads and 2 collars	SP-622-BC
	<b>2-0</b>	3	75	Blue	24	2 beads and 2 collars	SP-623-BC
	<b>2-0</b>	3	75	Blue	36		SP-623
<b>P-22</b> 1/2 Circle 13 mm 	<b>4-0</b>	1.5	45	Blue	12	Surgalloy™	SP-8656G


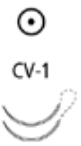

## TAPER POINT

Needles Specifications	Size USP	EP	Length cm	Color	Units Per box	Comments	Reorder code
<b>MV-70-3</b> 3/8 Circle 4 mm 	10-0	0.2	13	Blue	12	Surgalloy™	P-2790G
<b>MV-100-4</b> 3/8 Circle 5 mm 	10-0	0.2	13	Blue	12	Surgalloy™	P-2794G
<b>MV-100-4</b> 3/8 Circle 5 mm Double-arm 	10-0	0.2	5	Blue	12	Surgalloy™	P-2795G
<b>MV-135-5</b> 3/8 Circle 6 mm Double-arm 	8-0	0.4	45	Blue	36	Surgalloy™	VP-900-X
	8-0	0.4	5 × 45	Blue	12	Surgalloy™, LONG PACK	VP-901-MX
	8-0	0.4	2 × 45	Blue	12	Surgalloy™, LONG PACK	VP-900-MX
	8-0	0.4	60	Blue	36	Surgalloy™	VP-902-X
	8-0	0.4	5 × 60	Blue	12	Surgalloy™, LONG PACK	VP-903-MX
	8-0	0.4	2 × 60	Blue	12	Surgalloy™, LONG PACK	VP-902-MX
	8-0	0.4	75	Blue	12	Surgalloy™, LONG PACK	VP-973-X
	7-0	0.5	60	Blue	36	Surgalloy™	VP-904-X
	7-0	0.5	2 × 60	Blue	12	Surgalloy™, LONG PACK	VP-904-MX







## TAPER POINT

Needles Specifications	Size USP	EP	Length cm	Color	Units Per box	Comments	Reorder code
<b>MVF-135-5</b> 3/8 Circle 6 mm 	<b>8-0</b>	0.4	2 × 60	Blue	12	Surgalloy™, LONG PACK	VP-209MX
<b>MV-175-8</b> 3/8 Circle 8 mm Double-arm 	<b>8-0</b>	0.4	45	Blue	36	Surgalloy™	VP-744-X
	<b>8-0</b>	0.4	2 × 45	Blue	12	Surgalloy™, LONG PACK	VP-744-MX
	<b>8-0</b>	0.4	60	Blue	36	Surgalloy™	VP-745-X
	<b>8-0</b>	0.4	2 × 60	Blue	12	Surgalloy™, LONG PACK	VP-743-MX
	<b>8-0</b>	0.4	2 × 75	Blue	12	Surgalloy™, LONG PACK	VP-753-MX
	<b>7-0</b>	0.5	2 × 45	Blue	12	LONG PACK	VP-736-M
	<b>7-0</b>	0.5	45	Blue	36	Surgalloy™	VP-628
	<b>7-0</b>	0.5	5 × 45	Blue	12	Surgalloy™, LONG PACK	VP-734-MX
	<b>7-0</b>	0.5	60	Blue	36	Surgalloy™	VP-735-X
	<b>7-0</b>	0.5	2 × 60	Blue	12	Surgalloy™, LONG PACK	VP-737-MX
	<b>7-0</b>	0.5	5 × 60	Blue	12	Surgalloy™, LONG PACK	VP-735-MX
	<b>7-0</b>	0.5	75	Blue	36	Surgalloy™	VP-747-X
	<b>7-0</b>	0.5	2 × 75	Blue	12	Surgalloy™, LONG PACK	VP-747-MX
<b>MVF-175-8</b> 3/8 Circle 8 mm Double-arm 	<b>8-0</b>	0.4	2 × 60	Blue	12	Surgalloy™	VPF-738-MX
	<b>7-0</b>	0.5	60	Blue	36	Surgalloy™	VPF-735-X

## TAPER POINT







Needles Specifications	Size USP	EP	Length cm	Color	Units Per box	Comments	Reorder code	
<b>CV-351</b> 3/8 Circle 8 mm Double-arm		<b>8-0</b>	0.4	60	Blue	36	Cardioint™	VP-738-X
		<b>8-0</b>	0.4	2 × 60	Blue	12	Cardioint™, LONG PACK	VP-738-MX
		<b>7-0</b>	0.5	60	Blue	36	Cardioint™	VP-630-X
		<b>6-0</b>	0.7	60	Blue	36	Cardioint™	VP-631-X
<b>CV-1</b> 3/8 Circle 9 mm Double-arm		<b>7-0</b>	0.5	60	Blue	36	Optivis™	VP-72-X
		<b>7-0</b>	0.5	5 × 60	Blue	12	Optivis™, LONG PACK	VP-72-MX
		<b>7-0</b>	0.5	75	Blue	12	Optivis™, LONG PACK	VP-37-MX
		<b>7-0</b>	0.5	45	Blue	36	Surgalloy™, NUPACK™	VP-701-X
		<b>7-0</b>	0.5	5 × 45	Blue	12	Surgalloy™, LONG PACK	VP-701-MX
		<b>7-0</b>	0.5	60	Blue	36	Surgalloy™, NUPACK™	VP-702-X
		<b>7-0</b>	0.5	5 × 60	Blue	12	Surgalloy™, LONG PACK	VP-702-MX
		<b>7-0</b>	0.5	75	Blue	36	Surgalloy™, NUPACK™	VP-713-X
		<b>7-0</b>	0.5	5 × 75	Blue	12	Surgalloy™, LONG PACK	VP-713-MX
		<b>6-0</b>	0.7	45	Blue	36	Surgalloy™, NUPACK™	VP-746-X
		<b>6-0</b>	0.7	2 × 45	Blue	12	Surgalloy™, LONG PACK	VP-746-MX
		<b>6-0</b>	0.7	60	Blue	36	Surgalloy™, NUPACK™	VP-703-X
		<b>6-0</b>	0.7	60	Blue	12	Surgalloy™, LONG PACK	VP-703-MX
		<b>6-0</b>	0.7	75	Blue	12	Surgalloy™, LONG PACK	VP-709-MX
		<b>6-0</b>	0.7	75	Blue	36	Surgalloy™, NUPACK™	VP-709-X
		<b>6-0</b>	0.7	90	Blue	36	Surgalloy™, NUPACK™	VP-365-X
		<b>6-0</b>	0.7	60	Blue	36	Optivis™, NUPACK™	VP-73-X
		<b>6-0</b>	0.7	5 × 60	Blue	12	Optivis™, LONG PACK	VP-73-MX
<b>CVF-1</b> 3/8 Circle 9 mm Double-arm		<b>7-0</b>	0.5	2 × 60	Blue	12	Surgalloy™, LONG PACK	VP-102-MX
		<b>7-0</b>	0.5	60	Blue	36	Surgalloy™, NUPACK™	VPF-702-X
		<b>7-0</b>	0.5	75	Blue	36	Surgalloy™, NUPACK™	VPF-713-X
		<b>6-0</b>	0.7	60	Blue	36	Surgalloy™, NUPACK™	VPF-703-X
		<b>6-0</b>	0.7	75	Blue	36	Surgalloy™, NUPACK™	VPF-709-X
		<b>6-0</b>	0.7	5 × 60	Blue	12	Surgalloy™, LONG PACK	VPF-703-MX

## TAPER POINT







Needles Specifications	Size USP	EP	Length cm	Color	Units Per box	Comments	Reorder code
<b>MV-175-9</b> 3/8 Circle 10 mm Double-arm 	<b>8-0</b>	0.4	60	Blue	36	Surgalloy™	VP-763-X
	<b>7-0</b>	0.5	5 × 75	Blue	12	Surgalloy™, LONG PACK	VP-875-MX
	<b>7-0</b>	0.5	75	Blue	36	Surgalloy™	VP-875-X
	<b>7-0</b>	0.5	5 × 60	Blue	12	Surgalloy™, LONG PACK	VP-759-MX
<b>CV</b> 3/8 Circle 11 mm Double-arm 	<b>6-0</b>	0.7	5 × 60	Blue	36	Optivis™, NUPACK™	VP-74-X
	<b>6-0</b>	0.7	75	Blue	36	Surgalloy™, NUPACK	VP-776-X
<b>CVF</b> 3/8 Circle 11 mm Double-arm 	<b>6-0</b>	0.7	75	Blue	36	Surgalloy™, NUPACK	VPF-776-X
<b>CVL</b> Cardio curve 11 mm Double-arm 	<b>7-0</b>	0.5	2 × 60	Blue	12	Optivis™	VP-652-MX
<b>CV-11</b> 3/8 Circle 13 mm 	<b>6-0</b>	0.7	75	Blue	36	Surgalloy™	VP-889-X
	<b>5-0</b>	1	75	Blue	36	Surgalloy™	VP-890-X
<b>CV-11</b> 3/8 Circle 13 mm Double-arm 	<b>7-0</b>	0.5	75	Blue	36	Surgalloy™	VP-812-X
	<b>7-0</b>	0.5	60	Blue	36	Surgalloy™	VP-727-X
	<b>7-0</b>	0.5	45	Blue	36	Surgalloy™	VP-800-X
	<b>6-0</b>	0.7	90	Blue	36	Surgalloy™	VP-363X
	<b>6-0</b>	0.7	5 × 75	Blue	12	Surgalloy™	VP-706-MX
	<b>6-0</b>	0.7	75	Blue	36	Surgalloy™	VP-706-X
	<b>6-0</b>	0.7	5 × 45	Blue	12	Surgalloy™	VP-718-MX
	<b>6-0</b>	0.7	45	Blue	36	Surgalloy™	VP-718-X
	<b>6-0</b>	0.7	5 × 60	Blue	12	Surgalloy™	VP-726-MX
	<b>6-0</b>	0.7	60	Blue	36	Surgalloy™	VP-726-X
	<b>6-0</b>	0.7	5 × 75	Blue	12	Optivis™	VP-76-MX
	<b>6-0</b>	0.7	75	Blue	36	Optivis™	VP-76-X



## TAPER POINT

Needles Specifications	Size USP	EP	Length cm	Color	Units Per box	Comments	Reorder code	
<b>CV-11</b> 3/8 Circle 13 mm Double-arm		6-0	0.7	45	Blue	36	Optivis™	VP-78-X
		5-0	1	5 × 75	Blue	12	Surgalloy™	VP-705-MX
		5-0	1	75	Blue	36	Surgalloy™	VP-705-X
		5-0	1	90	Blue	36	Optivis™	VP-70-X
		5-0	1	45	Blue	36	Surgalloy™	VP-719-X
		5-0	1	90	Blue	36	Surgalloy™	VP-720-X
		5-0	1	5 × 60	Blue	12	Surgalloy™	VP-725-MX
		5-0	1	60	Blue	36	Surgalloy™	VP-725-X
		5-0	1	60	Blue	36	Optivis™	VP-75-X
		5-0	1	5 × 45	Blue	12	Surgalloy™ LONG PACK	VP-719-MX
		5-0	1	2 × 75	Blue	12	Surgalloy™ LONG PACK	VP-897-M
		5-0	1	5 × 90	Blue	12	Surgalloy™ LONG PACK	VP-720-MX
		4-0	1.5	75	Blue	36	Surgalloy™	VP-813-X
		<b>CVF-11</b> 3/8 Circle 13 mm Double-arm		6-0	0.7	60	Blue	36
6-0	0.7			5 × 60	Blue	12	Surgalloy™, LONG PACK	VP-126-MX
6-0	0.7			75	Blue	36	Surgalloy™, NUPACK	VPF-706-X
5-0	1			90	Blue	36	Surgalloy™	VPF-720-X
5-0	1			5 × 90	Blue	12	Surgalloy™, LONG PACK	VP-120-MX
<b>CV-301</b> 3/8 Circle 13 mm Double-arm		6-0	0.7	75	Blue	12		VP-86-X
<b>CV-15</b> 3/8 Circle 17 mm		4-0	1.5	75	Blue	36	Surgalloy™	VP-881-X
<b>CV-15</b> 3/8 Circle 17 mm Double-arm		5-0	1	75	Blue	36	Surgalloy™	VP-570X
		5-0	1	90	Blue	36	Surgalloy™, NUPACK	VP-580- X
		4-0	1.5	60	Blue	36	Surgalloy™	VP-584-X
		4-0	1.5	90	Blue	36	3 mm × 7 mm FIRM	VP-582-X
		4-0	1.5	90	Blue	36	Surgalloy™	VP-581-X
		4-0	1.5	120	Blue	36	Surgalloy™	VP-583-X
		3-0	2	90	Blue	36	Surgalloy™	VP-585-X
<b>CV-13</b> 3/8 Circle 22 mm Double-arm		4-0	1.5	75	Blue	36	Surgalloy™	VP-841-X
		4-0	1.5	90	Blue	36	Surgalloy™	VP-851-X
		3-0	2	90	Blue	36	Surgalloy™	VP-532-X
		2-0	3	90	Blue	36		VP-533






## TAPER POINT

Needles Specifications	Size	USP	EP	Length cm	Color	Units Per box	Comments	Reorder code
<b>CV-20</b> 1/2 Circle 10 mm Double-arm		<b>7-0</b>	0.5	60	Blue	36	Surgalloy™, NUPACK™	VP-772-X
		<b>6-0</b>	0.7	60	Blue	36	Surgalloy™, NUPACK™	VP-754-X
		<b>5-0</b>	1	60	Blue	36	Surgalloy™, NUPACK™	VP-755-X
<b>CV-337</b> 1/2 Circle 10 mm Double-arm		<b>6-0</b>	0.7	60	Blue	36	Cardiopoint™	VP-854-X
<b>CVF-21</b> 1/2 Circle 12 mm Double-arm		<b>4-0</b>	1.5	60	Blue	36	Surgalloy™	VPF-610-X
<b>CV-22</b> 1/2 Circle 13 mm Double-arm		<b>6-0</b>	0.7	45	Blue	36	Surgalloy™, NUPACK™	VP-714-X
		<b>6-0</b>	0.7	60	Blue	36	Surgalloy™, NUPACK™	VP-711-X
		<b>6-0</b>	0.7	5 × 75	Blue	12	Surgalloy™, LONG PACK	VP-733-MX
		<b>6-0</b>	0.7	75	Blue	36	Surgalloy™, NUPACK	VP-733-X
		<b>5-0</b>	1	45	Blue	36	Surgalloy™, NUPACK™	VP-715-X
		<b>5-0</b>	1	60	Blue	36	Surgalloy™, NUPACK™	VP-205-X
		<b>5-0</b>	1	5 × 75	Blue	12	Surgalloy™, LONG PACK	VP-710-MX
		<b>5-0</b>	1	75	Blue	36	Surgalloy™, NUPACK	VP-710-X
		<b>4-0</b>	1.5	60	Blue	36	Surgalloy™	VP-204-X
<b>CVF-22</b> 1/2 Circle 13 mm Double-arm		<b>6-0</b>	0.7	75	Blue	36	Surgalloy™, NUPACK	VPF-733-X
		<b>5-0</b>	1	75	Blue	36	Surgalloy™, NUPACK	VPF-710-X
		<b>5-0</b>	1	90	Blue	36	Surgalloy™, NUPACK	VPF-711-X
		<b>3-0</b>	2	90	Blue	12		VP-972X
<b>CV-23</b> 1/2 Circle 17 mm		<b>5-0</b>	1	75	Blue	36	Surgalloy™	VP-870-X
		<b>4-0</b>	1.5	75	Blue	36	Surgalloy™	VP-871-X
		<b>3-0</b>	2	75	Blue	36	Surgalloy™	VP-872-X

## TAPER POINT

Needles Specifications	Size USP	EP	Length cm	Color	Units Per box	Comments	Reorder code	
<b>CV-23</b> 1/2 Circle 17 mm Double-arm		6-0	0.7	90	Blue	36	Surgalloy™	VP-555X
		5-0	1	90	Blue	36	Optivis™	VP-660-X
		5-0	1	90	Blue	36	Surgalloy™	VP-556X
		5-0	1	90	Blue	36	Surgalloy™	VP-566X
		5-0	1	5 × 90	Blue	12	Surgalloy™	VP-556MX
		4-0	1.5	90	Blue	36	Surgalloy™	VPP-557-X
		4-0	1.5	5 × 90	Blue	12	Surgalloy™	VP-557MX
		4-0	1.5	45	Blue	36	Surgalloy™	VP-757-X
		4-0	1.5	90	Blue	36	Surgalloy™	VP-557X
		3-0	2	90	Blue	36	Surgalloy™	VP-558X
		2-0	3	90	Blue	36	455	VP-559-X
<b>CVF-23</b> 1/2 Circle 17 mm Double-arm		4-0	1.5	90	Blue	36	Surgalloy™	VPF-557-X
<b>CV-24</b> 1/2 Circle 20 mm Double-arm		4-0	1.5	90	Blue	36		VP-331X
<b>CV-25</b> 1/2 Circle 22 mm		4-0	1.5	75	Blue	36	Surgalloy™	VP-771-X
<b>CV-25</b> 1/2 Circle 22 mm Double-arm		4-0	1.5	90	Blue	36	Surgalloy™	VP-761-X
		4-0	1.5	5 × 90	Blue	12	Surgalloy™, LONG PACK	VP-761-MX
		3-0	2	75	Blue	36	Surgalloy™	VP-760-X
		3-0	2	90	Blue	36	Surgalloy™	VP-762-X

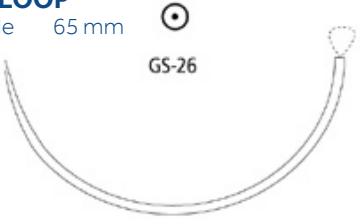
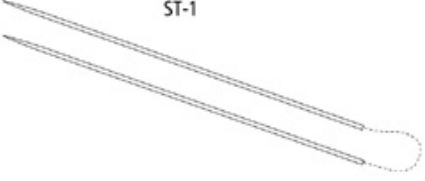
## TAPER POINT

Needles Specifications	Size		Length cm	Color	Units Per box	Comments	Reorder code
	USP	EP					
<b>V-20</b> 1/2 Circle 26 mm 	<b>4-0</b>	1.5	75	Blue	36	Surgalloy™	VP-831X
	<b>3-0</b>	2	75	Blue	36	Surgalloy™	VP-832X
	<b>2-0</b>	3	75	Blue	36	Surgalloy™	VP-833
	<b>0</b>	3.5	75	Blue	36	Surgalloy™	VP-834
<b>V-20</b> 1/2 Circle 26 mm Double-arm 	<b>5-0</b>	1	75	Blue	36	Surgalloy™, NUPACK	VP-510X
	<b>5-0</b>	1	90	Blue	36	Surgalloy™	VP-520X
	<b>4-0</b>	1.5	75	Blue	36	Surgalloy™	VP-511-X
	<b>4-0</b>	1.5	90	Blue	36	4.8 mm × 9.5 mm SOFT	VP-621-X
	<b>4-0</b>	1.5	90	Blue	36	Surgalloy™	VP-521-X
	<b>4-0</b>	1.5	90	Blue	36	3 mm × 7 mm FIRM	VP-421X
	<b>4-0</b>	1.5	120	Blue	36	Surgalloy™	VP-541X
	<b>3-0</b>	2	90	Blue	36	Surgalloy™	VP-522-X
	<b>3-0</b>	2	90	Blue	36	3 mm × 7 mm FIRM	VP-422X
	<b>3-0</b>	2	120	Blue	36	Surgalloy™	VP-542-X
	<b>2-0</b>	3	90	Blue	36	Surgalloy™	VP-523
	<b>2-0</b>	3	120	Blue	36	Surgalloy™	VP-543
	<b>VF-20</b> 1/2 Circle 26 mm 	<b>3-0</b>	2	45	Blue	36	Surgalloy™
<b>VF-20</b> 1/2 Circle 26 mm Double-arm 	<b>4-0</b>	1.5	90	Blue	36	Surgalloy™	VPF-521-X
	<b>3-0</b>	2	90	Blue	36	Surgalloy™	VPF-522-X
	<b>3-0</b>	2	135	Blue	36	Surgalloy™	VPF-354-X
	<b>2-0</b>	3	120	Blue	36	Surgalloy™	VPF-543
<b>GS-22</b> 1/2 Circle 27 mm 	<b>2-0</b>	3	5 × 45	Blue	12	D-Tach™	CP-16-M
	<b>2-0</b>	3	75	Blue	36		CP-411
	<b>0</b>	3.5	5 × 45	Blue	12	D-Tach™	CP-17-M
	<b>0</b>	3.5	75	Blue	36		CP-412



## TAPER POINT

Needles Specifications			Size		Length cm	Color	Units Per box	Comments	Reorder code
			USP	EP					
<b>V-30</b> 1/2 Circle    30 mm 		<b>2-0</b>	3	75	Blue	36		VP-945	
<b>GS-21</b> 1/2 Circle    37 mm 		<b>3-0</b>	2	75	Blue	36		CP-422	
		<b>2-0</b>	3	75	Blue	36		CP-423	
		<b>0</b>	3.5	5 × 45	Blue	12	D-Tach™	CP-10-M	
		<b>0</b>	3.5	75	Blue	36		CP-424	
		<b>1</b>	4	75	Blue	36		CP-425	
<b>V-26</b> 1/2 Circle    37 mm Double-arm 		<b>3-0</b>	2	90	Blue	36		VP-842-X	
		<b>2-0</b>	3	90	Blue	36		VP-843	
		<b>2-0</b>	3	120	Blue	36		VP-853	
<b>GS-24</b> 1/2 Circle    40 mm 		<b>2-0</b>	3	75	Blue	36		CP-433	
		<b>0</b>	3.5	75	Blue	36		CP-434	
		<b>1</b>	4	75	Blue	36		CP-435	
		<b>1</b>	4	100	Blue	36		CP-535	
<b>GS-25</b> 1/2 Circle    48 mm 		<b>0</b>	3.5	5 × 45	Blue	12	D-Tach™, LONG PACK	CP-30-M	
		<b>0</b>	3.5	75	Blue	36		CP-454	
		<b>1</b>	4	75	Blue	36		CP-455	
<b>GS-26</b> 1/2 Circle    65 mm 		<b>0</b>	3.5	3 × 60	Blue	12		3CP-444	
		<b>1</b>	4	150	Blue	24		CP-824	
		<b>2</b>	5	150	Blue	24		CP-825	






## TAPER POINT

Needles Specifications	Size		Length cm	Color	Units Per box	Comments	Reorder code
	USP	EP					
<b>GS-26 LOOP</b> 1/2 Circle 65 mm 	<b>1</b>	4	120	Blue	24	Loop	CP-824-L
	<b>1</b>	4	240	Blue	24	Loop	CP-834-L
<b>ST-1</b> Straight 51 mm Double-arm 	<b>2-0</b>	3	90	Blue	36		GP-282


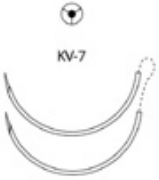
## HEAVY TAPER POINT

Needles Specifications	Size		Length cm	Color	Units Per box	Comments	Reorder code
	USP	EP					
<b>HGS-23</b> 1/2 Circle 22 mm 	<b>0</b>	3.5	5 × 45	Blue	12	D-Tach™, LONG PACK	CP-92-M
<b>HGS-22</b> 1/2 Circle 27 mm 	<b>0</b>	3.5	75	Blue	36		CP-414
	<b>1</b>	4	90	Blue	36		CP-415

## TAPERCUTTING

Needles Specifications	Size	USP	EP	Length cm	Color	Units Per box	Comments	Reorder code
<b>KV-1</b> 3/8 Circle 9 mm Double-arm		7-0	0.5	45	Blue	12	LONG PACK	VP-721-MX
		7-0	0.5	5 × 60	Blue	36	Surgalloy™	VP-704-X
		7-0	0.5	5 × 60	Blue	12	Surgalloy™, LONG PACK	VP-704-MX
		6-0	0.7	60	Blue	36	Surgalloy™, NUPACK	VP-712-X
		6-0	0.7	5 × 60	Blue	12	Surgalloy™, LONG PACK	VP-712-MX
<b>KVF-1</b> 3/8 Circle 9 mm Double-arm		7-0	0.5	60	Blue	36	Surgalloy™, NUPACK	VPF-704-X
		6-0	0.7	75	Blue	36	Surgalloy™, NUPACK	VPF-807-X
<b>YV-95</b> 3/8 Circle 9.5 mm Double-arm		7-0	0.5	60	Blue	36	Cardiopoint™	VP-911-X
<b>KV-11</b> 3/8 Circle 13 mm Double-arm		6-0	0.7	60	Blue	36	Surgalloy™, NUPACK	VP-707-X
		6-0	0.7	75	Blue	36	Surgalloy™, NUPACK	VP-728-X
		6-0	0.7	5 × 75	Blue	12	Surgalloy™, LONG PACK	VP-728-MX
		5-0	1	75	Blue	36	Surgalloy™, NUPACK	VP-708-X
		5-0	1	90	Blue	36	Surgalloy™, NUPACK	VP-751-X
<b>KVF-11</b> 3/8 Circle 13 mm Double-arm		6-0	0.7	75	Blue	36	Surgalloy™, NUPACK	VPF-728-X

## TAPERCUTTING

Needles Specifications	Size		Length cm	Color	Units Per box	Comments	Reorder code
	USP	EP					
<b>KV-15</b> 3/8 Circle 17 mm 	<b>5-0</b>	1	75	Blue	36	Surgalloy™, NUPACK	VP-941-X
	<b>4-0</b>	1.5	75	Blue	36	Surgalloy™	VP-942-X
<b>KV-15</b> 3/8 Circle 17 mm Double-arm 	<b>5-0</b>	1	75	Blue	36	Surgalloy™, NUPACK	VP-914-X
	<b>5-0</b>	1	90	Blue	36	Surgalloy™	VP-924-X
	<b>4-0</b>	1.5	90	Blue	36	Surgalloy™	VP-925-X
	<b>3-0</b>	2	90	Blue	36	Surgalloy™	VP-926-X
<b>KV-26</b> 3/8 Circle 26 mm Double-arm 	<b>4-0</b>	1.5	120	Blue	36	Surgalloy™	VP-965
	<b>3-0</b>	1.5	120	Blue	36		VP-956
<b>KV-5</b> 1/2 Circle 17 mm Double-arm 	<b>5-0</b>	1	90	Blue	36	Surgalloy™, NUPACK	VP-934-X
	<b>4-0</b>	1.5	90	Blue	36	Surgalloy™	VP-935-X
	<b>4-0</b>	1.5	90	Blue	36	3 mm × 7 mm FIRM	VP-435-X
	<b>3-0</b>	2	90	Blue	36	Surgalloy™	VP-936-X
	<b>2-0</b>	3	90	Blue	36	Surgalloy™	VP-937
<b>KV-7</b> 1/2 Circle 26 mm Double-arm 	<b>4-0</b>	1.5	90	Blue	36	Surgalloy™	VP-975-X
	<b>3-0</b>	2	90	Blue	36	Surgalloy™	VP-976-X
	<b>3-0</b>	2	90	Blue	36	3 mm × 7 mm FIRM	VP-986-X
	<b>2-0</b>	3	90	Blue	36	Surgalloy™	VP-977
	<b>0</b>	3.5	90	Blue	36	Surgalloy™	VP-978
<b>KV-34</b> 1/2 Circle 37 mm 	<b>0</b>	3.5	75	Blue	36		CP-844
	<b>1</b>	4	75	Blue	36		CP-845



## BLUNT POINT

Needles Specifications	Size		Length cm	Color	Units Per box	Comments	Reorder code
	USP	EP					
<b>BGS-21</b> 1/2 Circle 37 mm 	<b>0</b>	3.5	75	Blue	36	ProtectPoint™	CP-957
<b>BGS-25</b> 1/2 Circle 48 mm 	<b>1</b>	4	100	Blue	36	ProtectPoint™	CP-836

## SPATULA

Needles Specifications	Size		Length cm	Color	Units Per box	Comments	Reorder code
	USP	EP					
<b>SE-140-8</b> 3/8 Circle 6.3 mm Double-arm Curvature 140° 	<b>10-0</b>	0.2	30	Blue	12	Surgalloy™	P-2771K
<b>SE-160-6</b> 1/2 Circle 5.5 mm Double-arm Curvature 160° 	<b>10-0</b>	0.2	30	Blue	12	Surgalloy™	P-2756K

# TI-CRON™ / SURGIDAC™

Coated and Uncoated Braided and Monofilament Polyester

A brand you have trusted for over 35 years.



Please refer to the individual package IFUs for complete instructions, indications, contraindications, warnings and precautions.

<b>Description</b>	Ti-Cron™ and Surgidac™ polyester sutures are non absorbable, sterile, surgical sutures composed of Polyethylene Terephthalate. They are prepared from fibers of high molecular weight, long-chain, linear polyesters having recurrent aromatic rings as an integral component. Ti-Cron™ braided sutures are available coated uniformly with silicone to increase surface lubricity, thereby enhancing handling characteristics, ease of passage through tissue, and knot run-down properties. Surgidac™ monofilament sutures are uncoated and the braided sutures are available uncoated or coated uniformly with polybutylene adipate to increase surface lubricity, thereby enhancing handling characteristics, ease of passage through tissue, and knot run-down properties.
<b>Indications</b>	Ti-Cron™ and Surgidac™ polyester sutures are indicated for use in general soft tissue approximation and/or ligation, including use in cardiovascular, ophthalmic and neurological surgery.
<b>Actions</b>	Ti-Cron™ and Surgidac™ polyester sutures elicit a minimal acute inflammatory reaction in tissues, followed by gradual encapsulation of the suture by fibrous connective tissue.
<b>Tensile Strength</b>	Ti-Cron™ and Surgidac™ polyester sutures are not absorbed, nor is any significant change in tensile strength retention known to occur in vivo.
<b>Sterilization</b>	Ti-Cron™ and Surgidac™ polyester sutures in sizes 6-0 or smaller or pledgeted sutures are sterilized by ethylene oxide. Sizes 5-0 or greater are sterilized by gamma irradiation.
<b>Packaging</b>	Ti-Cron™ polyester sutures are available undyed (white) or dyed with D&C Blue No. 6 and braided, coated, in sizes 5 (7 Metric), size 2 (5 metric), size 1 (4 metric), 0 size (3.5 metric) through 5-0 (1 Metric). SURGIDAC™ polyester sutures are available undyed (white), braided and uncoated in USP and EP sizes 5-0 (1 metric). The sutures are supplied sterile, in pre-cut lengths, affixed to various needle types using both permanent and removable needle attachment techniques. Ti-Cron™ polyester sutures are also available with PTFE (polytetrafluorethylene) pledgets for use as a pad between the suture area and the tissue surface to increase the load bearing area. Please refer to the individual package IFUs for complete instructions, indications, contraindications, warnings and precautions.






## Ti-Cron™ / Surgidac™

Ti-Cron™ / Surgidac™ is a suture brand you have trusted for over 35 years with a proven history for performance


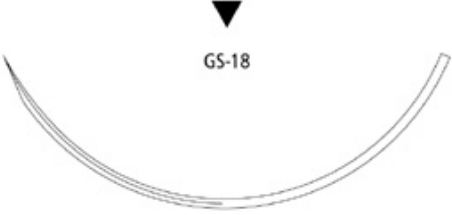
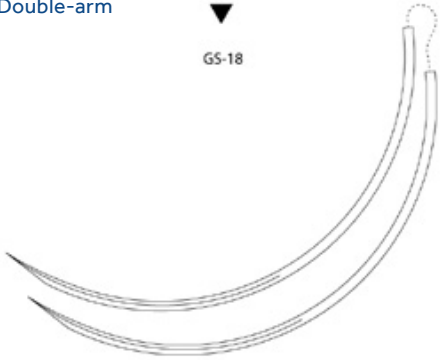



PRE-CUT							
Needles Specifications	Size		Length cm	Color	Units Per box	Comments	Reorder code
	USP	EP					
	5-0	1	12×45	Blue	36	Pre-cut	88863008-21
	4-0	1.5	12×45	Blue	36	Pre-cut	88863008-31
	4-0	1.5	12×75	Blue	24	Pre-cut	88863003-32
	3-0	2	12×45	Blue	36	Pre-cut	88863008-41
	3-0	2	6×45	Blue	36	Pre-cut	88863001-41
	3-0	2	12×75	Blue	24	Pre-cut	88863003-42
	2-0	3	5×45	White	12	Pre-cut	88863011-51
	2-0	3	12×45	Blue	36	Pre-cut	88863008-51
	2-0	3	12×75	Blue	24	Pre-cut	88863003-52
	0	3.5	6×45	Blue	24	Pre-cut	88863001-62
	0	3.5	6×75	Blue	24	Pre-cut	88863002-62
	1	4	6×45	Blue	24	Pre-cut	88863001-72

REEL							
Needles Specifications	Size		Length cm	Color	Units Per box	Comments	Reorder code
	USP	EP					
	3-0	2	250	Blue	24	Reel	88863380-42
	2-0	3	250	Blue	24	Reel	88863380-52
	0	3.5	250	Blue	24	Reel	88863380-62
	1	4	250	Blue	24	Reel	88863380-72
	2	5	250	Blue	24	Reel	88863380-82
	4	6	150	Blue	24	Reel	88863381-09
	5	7	150	Blue	24	Reel	88863381-89

## REVERSE CUTTING

Needles Specifications			Size		Length cm	Color	Units Per box	Comments	Reorder code
			USP	EP					
<b>C-12</b> 3/8 Circle 16 mm			5-0	1	45	Blue	36		88863045-21
			4-0	1.5	75	Blue	36		88863046-31
<b>C-13</b> 3/8 Circle 19 mm			4-0	1.5	45	Blue	36		88863047-31
			3-0	2	75	Blue	36		88863047-41
			3-0	2	45	Blue	36		88863055-41
			2-0	3	75	Blue	36		88863048-51
<b>C-14</b> 3/8 Circle 24 mm			3-0	2	45	Blue	36		88863069-41
			3-0	2	75	Blue	36		88863050-41
			2-0	3	75	Blue	36		88863050-51
			2-0	3	90	Blue	36		88863065-51
			0	3.5	75	Blue	36		88863050-61
			0	3.5	90	Blue	36		88863065-61
<b>C-15</b> 3/8 Circle 26 mm			2-0	3	75	Blue	36		88863070-51
<b>C-16</b> 3/8 Circle 30 mm			2-0	3	75	Blue	36		88863023-51
			0	3.5	75	Blue	36		88863023-61
			0	3.5	90	Blue	36		88863012-61
			1	4	75	Blue	36		88863023-71



## REVERSE CUTTING

Needles Specifications	Size		Length cm	Color	Units Per box	Comments	Reorder code
	USP	EP					
<b>C-17</b> 3/8 Circle    39 mm 	<b>2-0</b>	3	75	Blue	36		88863054-51
	<b>0</b>	3.5	75	Blue	36		88863054-61
	<b>1</b>	4	75	Blue	36		88863054-71
	<b>2</b>	5	75	Blue	36		88863054-81
<b>GS-18</b> 3/8 Circle    77 mm 	<b>5</b>	7	45	Blue	12		88863056-89
<b>GS-18</b> 3/8 Circle    77 mm Double-arm 	<b>5</b>	7	75	Blue	24		88863128-79
	<b>2</b>	5	75	Blue	24		88863128-82
	<b>2</b>	5	75	Green	24		CD-496
<b>HOS-10</b> 1/2 Circle    26 mm 	<b>2-0</b>	3	75	Blue	36		88863395-51
	<b>0</b>	3.5	75	Blue	36		88863395-61
	<b>1</b>	4	75	Blue	36		88863395-71
	<b>2</b>	5	75	Blue	36		88863395-81
<b>GS-11</b> 1/2 Circle    37 mm 	<b>2-0</b>	3	75	Blue	36		88863309-51
	<b>0</b>	3.5	75	Blue	36		88863309-61
	<b>1</b>	4	75	Blue	36		88863309-71
	<b>4</b>	6	75	Blue	36		88863154-09
<b>HOS-11</b> 1/2 Circle    37 mm 	<b>1</b>	4	75	Blue	36		88863028-71
	<b>2</b>	5	75	Blue	36		88863028-81

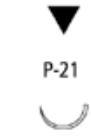
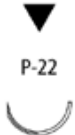

## REVERSE CUTTING

Needles Specifications	Size		Length cm	Color	Units Per box	Comments	Reorder code
	USP	EP					
<b>GS-12</b> 1/2 Circle    40 mm 	2	5	75	Blue	36		88863155-81
<b>HOS-12</b> 1/2 Circle    40 mm 	1	4	75	Blue	36		88863026-71
	2	5	75	Blue	36		88863026-81
<b>HOS-14</b> 1/2 Circle    55 mm 	5	7	75	Blue	36		88863027-79



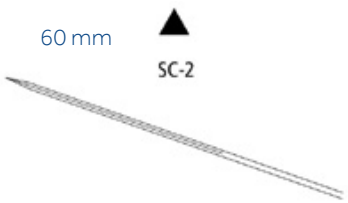
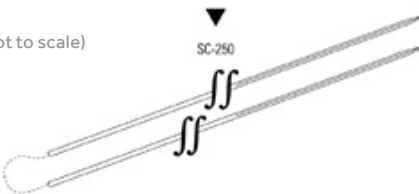
## PREMIUM REVERSE CUTTING

Needles Specifications	Size		Length cm	Color	Units Per box	Comments	Reorder code
	USP	EP					
<b>P-13</b> 3/8 Circle    13 mm 	5-0	1	45	Blue	12	Surgalloy™	SCD-3058G
	5-0	1	45	White	12	Surgalloy™	SCD-2864G
	4-0	1.5	45	Blue	12	Surgalloy™	SCD-3059G
	4-0	1.5	45	White	12	Surgalloy™	88862862-34
<b>P-12</b> 3/8 Circle    19 mm 	4-0	1.5	45	Blue	12	Surgalloy™	SCD-3049G

## PREMIUM REVERSE CUTTING


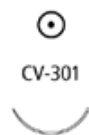






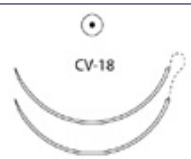
Needles Specifications	Size		Length cm	Color	Units Per box	Comments	Reorder code
	USP	EP					
<b>P-21</b> 1/2 Circle 9 mm 	<b>4-0</b>	1.5	45	White	12		88862854-33
<b>P-22</b> 1/2 Circle 13 mm 	<b>3-0</b>	2	45	Blue	12	Surgalloy™	88863109-43
<b>P-24</b> 1/2 Circle 16 mm 	<b>2-0</b>	3	45	Blue	12	Surgalloy™	88863059-53

## CONVENTIONAL CUTTING





Needles Specifications	Size		Length cm	Color	Units Per box	Comments	Reorder code
	USP	EP					
<b>SCC</b> 1/2 Circle 48 mm 	<b>5</b>	7	5×75	Blue	6		88863127-79
<b>SC-1</b> Straight 51 mm 	<b>2-0</b>	3	75	Blue	36		88863163-51
	<b>0</b>	3.5	75	Blue	36		88863163-61
	<b>1</b>	4	75	Blue	36		88863163-71
<b>SC-2</b> Straight 60 mm 	<b>3-0</b>	2	75	Blue	36		88863062-41
	<b>2-0</b>	3	75	Blue	36		88863062-51
	<b>2-0</b>	3	90	Blue	36		88863165-51
	<b>0</b>	3.5	75	Blue	36		88863062-61
<b>SC-250</b> Straight 254 mm (Not to scale) 	<b>2-0</b>	3	75	Blue	12		88862947-53





## TAPER POINT

Needles Specifications	Size USP	EP	Length cm	Color	Units Per box	Comments	Reorder code
<b>CV-310</b> 3/8 Circle 10 mm Double-arm	7-0	0.5	45	Blue	36	Cardiopoint™	88863288-01
							
<b>CV-301</b> 3/8 Circle 13 mm	4-0	1.5	75	Blue	36	Cardiopoint™	88863106-31
							
<b>CV-301</b> 3/8 Circle 13 mm Double-arm	6-0	0.7	75	Blue	36	Cardiopoint™	88863280-11
							
	5-0	1	75	Blue	36	3mm x 3mm FIRM, Cardiopoint™	88863314-21
	5-0	1	75	Blue	36	Cardiopoint™	88863280-21
	5-0	1	75	White	36	Cardiopoint™	88862801-21
	4-0	1.5	75	Blue	36	3 mm x 3 mm Firm, Cardiopoint™	88863314-31
	4-0	1.5	75	Blue	36	Cardiopoint™	88863280-31
3-0	2	75	Blue	36	Cardiopoint™	88863280-41	
							
<b>CV-15</b> 3/8 Circle 17 mm Double-arm	2-0	3	90	White	6	Surgalloy™	88862636-56
							
<b>CV-327</b> 3/8 Circle 18 mm Double-arm	5-0	1	75	Blue	36	Cardiopoint™	88862775-21
							
	4-0	1.5	75	Blue	36	Cardiopoint™	88862775-31
3-0	2	75	Blue	36	Cardiopoint™	88862775-41	
							
<b>CV-13</b> 3/8 Circle 22 mm Double-arm	3-0	2	90	Blue	36	Surgalloy™	88862627-41
							
<b>CV-18</b> 3/8 Circle 26 mm Double-arm	2-0	3	10x90	5 Blue and 5 White	6		88863230-56
							


## TAPER POINT

Needles Specifications	Size USP	EP	Length cm	Color	Units Per box	Comments	Reorder code	
<b>CV-337</b> 1/2 Circle 10 mm Double-arm		5-0	1	75	Blue	36	Cardiopoint™	88863298-21
		5-0	1	90	Blue	36	Pledgets 3 mm × 3 mm FIRM, Cardiopoint™	88863297-21
		5-0	1	90	Blue	36	Cardiopoint™	88863139-21
<b>CV-330</b> 1/2 Circle 13 mm Double-arm		5-0	1	60	Blue	36	Cardiopoint™	88863226-21
		4-0	1.5	60	Blue	36	Cardiopoint™	88863226-31
		4-0	1.5	90	Blue	36	Cardiopoint™	88863249-31
		3-0	2	45	Blue	36	Pledgets 3 mm × 3 mm FIRM, Cardiopoint™	88863295-41
		3-0	2	60	Blue	36	Cardiopoint™	88863226-41
<b>CV-331</b> 1/2 Circle 16 mm		5-0	1	75	Blue	36	Cardiopoint™	88863186-21
		4-0	1.5	75	Blue	36	Cardiopoint™	88863186-31
		3-0	2	75	Blue	36	Cardiopoint™	88863186-41
		2-0	3	75	Blue	36	Cardiopoint™	88863186-51
<b>CV-331</b> 1/2 Circle 16 mm Double-arm		5-0	1	90	Blue	36	Cardiopoint™	88863271-21
		4-0	1.5	75	Blue	36	Cardiopoint™	88863087-31
		4-0	1.5	90	Blue	36	Cardiopoint™	88863271-31
		3-0	2	10×75	Blue	6	Cardiopoint™	88863084-46
		3-0	2	75	Blue	36	Cardiopoint™	88863087-41
		3-0	2	90	Blue	36	Cardiopoint™	88863271-41
		2-0	3	8×75	4 Blue and 4 White	6	3 mm × 3 mm FIRM, Cardiopoint™	88863086-56
		2-0	3	8×75	4 Blue and 4 White	6	Cardiopoint™	88863083-56
		2-0	3	10×75	Blue	6	Cardiopoint™	88863084-56
		2-0	3	8×75	4 Blue and 4 White	6	Pledgets 3 mm × 7 mm FIRM, Cardiopoint™	88863091-56
		2-0	3	75	Blue	36	Cardiopoint™	88863087-51
		2-0	3	90	5 Blue and 5 White	6	Pledgets 3 mm × 3 mm FIRM, Cardiopoint™	88863273-56
		2-0	3	10×90	5 Blue and 5 White	6	Pledgets 3mm × 7 mm FIRM, Cardiopoint™	88863269-56

## TAPER POINT

Needles Specifications	Size		Length cm	Color	Units Per box	Comments	Reorder code
	USP	EP					
<b>CV-24</b> 1/2 Circle 20 mm 	<b>3-0</b>	2	5×45	Blue	24	D-Tach™	88863267-42
<b>CV-316</b> 1/2 Circle 20 mm Double-arm 	<b>4-0</b>	1.5	75	Blue	36	Pledgets 3 mm × 7 mm FIRM, Cardiopoint™	88863348-31
	<b>4-0</b>	1.5	90	Blue	36	Cardiopoint™	88863369-31
	<b>3-0</b>	2	75	Blue	36	Cardiopoint™	88863227-41
	<b>3-0</b>	2	90	Blue	36	Cardiopoint™	88863369-41
	<b>2-0</b>	3	10×75	5 Blue and 5 White	6	Cardiopoint™	88863385-56
	<b>2-0</b>	3	75	Blue	36	Cardiopoint™	88863227-51
	<b>2-0</b>	3	75	Blue	36	3 mm × 7 mm FIRM, Cardiopoint™	88863348-51
	<b>2-0</b>	3	10×75	5 Blue and 5 White	6	3 mm × 7 mm FIRM, Cardiopoint™	88863323-56
	<b>2-0</b>	3	75	White	36	3 mm × 7 mm FIRM, Cardiopoint™	88862978-51
	<b>2-0</b>	3	10×75	5 Blue and 5 White	6	Pledgets 3 mm × 3 mm FIRM, Cardiopoint™	88863277-56
	<b>2-0</b>	3	10×75	Blue	6	Cardiopoint™	88863344-56
	<b>2-0</b>	3	90	White	36	Cardiopoint™	88862839-51
	<b>2-0</b>	3	90	Blue	36	Cardiopoint™	88863369-51

## TAPER POINT

Needles Specifications	Size		Length cm	Color	Units Per box	Comments	Reorder code
	USP	EP					
<b>GS-23</b> 1/2 Circle 22 mm 	<b>0</b>	3.5	5×45	Blue	24	D-Tach™	88863257-62
<b>CV-25</b> 1/2 Circle 22 mm Double-arm 	<b>2-0</b>	3	90	Blue	36	Surgalloy™	88863370-51





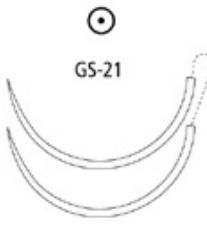
## TAPER POINT

Needles Specifications	Size		Length cm	Color	Units Per box	Comments	Reorder code
	USP	EP					
<b>CV-305</b> 1/2 Circle 25 mm 	4-0	1.5	75	Blue	36	Cardiopoint™	88863185-31
	3-0	2	75	Blue	36	Cardiopoint™	88863185-41
	2-0	3	10×75	5 Blue and 5 White	6	Cardiopoint™	88863353-56
	2-0	3	5×75	Blue	24	Cardiopoint™	88863033-52
	0	3.5	75	Blue	36	Cardiopoint™	88863185-61
	1	4	75	Blue	36	Cardiopoint™	88863185-71
	2-0	3	75	Blue	36	Cardiopoint™	88863185-51
	<b>CV-305</b> 1/2 Circle 25 mm Double-arm 	4-0	1.5	90	Blue	36	Cardiopoint™
3-0		2	90	Blue	36	Pledgets 3 mm × 7 mm FIRM, Cardiopoint™	88863336-41
3-0		2	90	Blue	36	Cardiopoint™	88863256-41
3-0		2	105	Blue	36	Cardiopoint™	88863035-41
2-0		3	10×75	5 Blue and 5 White	6	3 mm × 7 mm FIRM, Cardiopoint™	88863324-56
2-0		3	8×75	4 Blue AND 4 White	6	Pledgets 3 mm × 7 mm SOFT, Cardiopoint™	VCD-346
2-0		3	10×75	5 Blue and 5 White	6	Cardiopoint™	88863356-56
2-0		3	5×75	Blue	6	Cardiopoint™	88863040-56
2-0		3	10×75	5 Blue and 5 White	6	Pledgets 4 mm × 8 mm FIRM, Cardiopoint™	88863037-56
2-0		3	90	Blue	36	3 mm × 7 mm FIRM, Cardiopoint™	88863336-51
2-0		3	90	Blue	36	Cardiopoint™	88863256-51
2-0		3	10×90	5 Blue and 5 White	6	Cardiopoint™	88862981-56
2-0		3	90	White	36	Cardiopoint™	88862914-51
2-0		3	90	White	36	3 mm × 7 mm FIRM, Cardiopoint™	88862953-51

## TAPER POINT

Needles Specifications	Size		Length cm	Color	Units Per box	Comments	Reorder code
	USP	EP					
<b>CV-305</b> 1/2 Circle 25 mm Double-arm 	<b>2-0</b>	3	8×90	4 Blue and 4 White	6	Pledgets 3 mm × 3 mm FIRM, Cardiopoint™	88863307-56
	<b>2-0</b>	3	6×90	3 Blue and 3 White	6	Pledgets 3 mm × 7 mm FIRM, Cardiopoint™	88863306-56
	<b>2-0</b>	3	10×90	Blue	6	Cardiopoint™	88863220-56
	<b>2-0</b>	3	105	Blue	36	Cardiopoint™	88863035-51
	<b>0</b>	3.5	90	Blue	36	Cardiopoint™	88863256-61
	<b>0</b>	3.5	90	White	36	Cardiopoint™	88862914-61
<b>V-20</b> 1/2 Circle 26 mm 	<b>4-0</b>	1.5	5×45	Blue	24	D-Tach™	88863258-32
	<b>3-0</b>	2	5×45	Blue	24	D-Tach™	88863258-42
	<b>2-0</b>	3	5×45	Blue	24	D-Tach™	88863258-52
	<b>2-0</b>	3	120	Blue	36	Endoscopy	88863129-51
	<b>0</b>	3.5	120	Blue	36	Endoscopy	88863015-61
<b>V-20</b> 1/2 Circle 26 mm Double-arm 	<b>2-0</b>	3	90	Blue	36	Surgalloy™	88863251-51
<b>GS-22</b> 1/2 Circle 27 mm 	<b>2-0</b>	3	75	Blue	36		88863093-51
	<b>0</b>	3.5	5×45	Blue	24	D-Tach™	88863259-62
	<b>0</b>	3.5	75	Blue	36		88863093-61

## TAPER POINT





Needles Specifications		Size USP	EP	Length cm	Color	Units Per box	Comments	Reorder code
<b>V-30</b> 1/2 Circle 30 mm		<b>3-0</b>	2	75	Blue	36		88863205-41
		<b>2-0</b>	3	75	Blue	36		88863205-51
		<b>0</b>	3.5	75	Blue	36		88863205-61
<b>CV-300</b> 1/2 Circle 35 mm		<b>0</b>	3.5	75	Blue	36	Cardiopoint™	88863082-61
		<b>2-0</b>	3	75	Blue	36	Cardiopoint™	88863082-51
<b>CV-300</b> 1/2 Circle 35 mm Double-arm		<b>2-0</b>	3	75	Blue	36	Cardiopoint™	88863274-51
		<b>2-0</b>	3	90	Blue	36	Cardiopoint™	88863393-51
		<b>0</b>	3.5	75	Blue	36	Cardiopoint™	88863274-61
		<b>0</b>	3.5	90	Blue	36	Cardiopoint™	88863393-61
<b>GS-21</b> 1/2 Circle 37 mm		<b>2-0</b>	3	5×45	Blue	24	D-Tach™	88863260-52
		<b>2-0</b>	3	75	Blue	36		88863090-51
		<b>0</b>	3.5	45	Blue	24	D-Tach™	88863260-62
		<b>0</b>	3.5	75	Blue	36		88863090-61
		<b>1</b>	4	5×45	Blue	24	D-Tach™	88863260-72
		<b>1</b>	4	75	Blue	36		88863090-71
<b>GS-21</b> 1/2 Circle 37 mm Double-arm		<b>2</b>	5	90	Blue	36		88863146-81
		<b>2</b>	5	75	Blue	36		88863090-81

## TAPER POINT


Needles Specifications	Size		Length cm	Color	Units Per box	Comments	Reorder code
	USP	EP					
<b>GS-25</b> 1/2 Circle 48 mm 	0	3.5	5×45	Blue	24	D-Tach™	88863261-62
	0	3.5	75	Blue	36		88863097-61
	1	4	5×45	Blue	24	D-Tach™	88863261-72
	1	4	75	Blue	36		88863097-71
	2	5	5×45	Blue	24	D-Tach™	88863261-82
<b>GS-27</b> 1/2 Circle 76 mm Double-arm 	0	3.5	75	Blue	36		88863190-61
<b>GU-46</b> 5/8 Circle 27 mm 	2-0	3	75	Blue	36		88862768-51
<b>SK</b> SKI 23 mm 	2-0	3	120	Blue	36		88863017-51
	0	3.5	120	Blue	36		88863017-61
<b>TJ-32</b> 1/2 Circle 32 mm 	1	4	75	Blue	36		88863103-71





## HEAVY TAPER POINT

Needles Specifications	Size USP	EP	Length cm	Color	Units Per box	Comments	Reorder code
<b>HGS-23</b> 1/2 Circle    22 mm 	<b>2-0</b>	3	5×45	Blue	24	D-Tach™	88863374-52
	<b>0</b>	3.5	5×45	Blue	24	D-Tach™	88863374-62
	<b>1</b>	4	75	Blue	36		88863375-71
<b>HGU-46</b> 5/8 Circle    27 mm 	<b>0</b>	3.5	5×45	Blue	24	D-Tach™	88863386-62
<b>HGS-22</b> 1/2 Circle    27 mm 	<b>2-0</b>	3	5×45	Blue	24	D-Tach™	88863379-52
	<b>1</b>	4	5×45	Blue	24	D-Tach™	88863372-72
	<b>2</b>	5	5×45	Blue	24		88863372-82
<b>HGS-21</b> 1/2 Circle    37 mm 	<b>1</b>	4	75	Blue	36		88863113-71
	<b>2</b>	5	75	Blue	36		88863113-81
	<b>5</b>	7	75	Blue	36		88863111-79


## HEAVY TAPER POINT

Needles Specifications	Size		Length cm	Color	Units Per box	Comments	Reorder code
	USP	EP					
<b>HGS-24</b> 1/2 Circle    40 mm 	<b>2</b>	5	75	Blue	36		88863376-81

## PENETRATING TAPER

Needles Specifications	Size		Length cm	Color	Units Per box	Comments	Reorder code
	USP	EP					
<b>YE-7</b> 3/8 Circle    18 mm 	<b>4-0</b>	1.5	75	Blue	36	Cardiopoint™	88863107-31
<b>YE-7</b> 3/8 Circle    18 mm <b>Double-arm</b> 	<b>4-0</b>	1.5	75	Blue	36	Cardiopoint™	88863182-31
	<b>3-0</b>	2	75	Blue	36	Cardiopoint™	88863182-41
	<b>2-0</b>	3	10×75	5 Blue and 5 White	6	Pledgets 3 mm × 7 mm FIRM, Cardiopoint™	88863365-56
	<b>2-0</b>	3	75	Blue	36	Pledgets 3 mm × 7 mm FIRM, Cardiopoint™	88863233-51
	<b>2-0</b>	3	10×75	5 Blue and 5 White	6	Cardiopoint™	88863211-56
	<b>2-0</b>	3	75	Blue	36	Cardiopoint™	88863182-51
	<b>2-0</b>	3	5×90	5 Blue and 5 White	6	Cardiopoint™	88863052-56
	<b>2-0</b>	3	90	5 Blue and 5 White	6	Pledgets 3 mm × 7 mm FIRM, Cardiopoint™	88862658-56
	<b>2-0</b>	3	90	Blue	36	Cardiopoint™	88862631-51


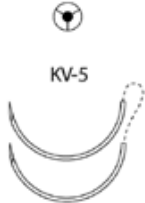

## PENETRATING TAPER

Needles Specifications	Size		Length cm	Color	Units Per box	Comments	Reorder code	
	USP	EP						
<b>Y-31</b> 1/2 Circle 16 mm Double-arm 	4-0	1.5	8×75	4 Blue and 4 White	6	Pledgets 3 mm × 3 mm FIRM, Cardioint™	88863218-36	
	4-0	1.5	90	Blue	36	Cardioint™	88863160-31	
	3-0	2	75	White	36	Cardioint™	88862880-41	
	3-0	2	90	Blue	36	Cardioint™	88863160-41	
	2-0	3	10×75	5 Blue and 5 White	6	3 mm v 7 mm FIRM, Cardiointv	88863221-56	
	2-0	3	8×75	4 Blue and 4 White	6	3 mm × 7 mm FIRM, Cardiointv	88863219-56	
	2-0	3	8×75	4 Blue and 4 White	6	3 mm × 3 mm FIRM, Cardiointv	88863218-56	
	2-0	3	75	Blue	36	Cardioint™	88863125-51	
	2-0	3	75	White	36	Cardiointv	88862880-51	
	2-0	3	10×75	5 Blue and 5 White	6	Cardiointv	88862877-56	
	2-0	3	90	White	36	3 mm × 3 mm FIRM, Cardiointv	88862811-51	
	2-0	3	90	White	36	3 mm × 7 mm FIRM, Cardioint™	88862808-51	
	2-0	3	90	Blue	36	3 mm × 3 mm FIRM, Cardioint™	88863191-51	
	2-0	3	10×90	5 Blue and 5 White	6	Pledgets 3 mm × 3 mm FIRM, Cardioint™	88863255-56	
	2-0	3	90	Blue	36	Pledgets 3 mm × 7 mm FIRM, Cardioint™	88863204-51	
	2-0	3	90	Blue	36	Cardioint™	88863160-51	
	2-0	3	10×90	5 Blue and 5 White	6	Cardioint™	88862851-56	
	<b>Y-16</b> 1/2 Circle 20 mm Double-arm 	2-0	3	6×75	3 Blue and 3 White	6	Pledgets 3 mm × 7 mm FIRM, Cardioint™	88863303-56
		2-0	3	10×75	5 Blue and 5 White	6	Cardioint™	88863241-56
2-0		3	90	Blue	36	Cardioint™	88863157-51	



## PENETRATING TAPER

Needles Specifications	Size		Length cm	Color	Units Per box	Comments	Reorder code
	USP	EP					
<b>Y-5</b> 1/2 Circle 26 mm 	<b>3-0</b>	2	75	Blue	36	Cardiopoint™	88863088-41
	<b>2-0</b>	3	75	Blue	36	Cardiopoint™	88863088-51
<b>Y-5</b> 1/2 Circle 26 mm Double-arm 	<b>4-0</b>	1.5	90	Blue	36	Cardiopoint™	88863167-31
	<b>3-0</b>	2	75	Blue	36	Cardiopoint™	88863088-41
	<b>3-0</b>	2	90	Blue	36	Cardiopoint™	88863167-41
	<b>2-0</b>	3	10×75	5 Blue and 5 White	6	Cardiopoint™	88863222-56
	<b>2-0</b>	3	8×75	4 Blue and 4 White	6	Pledgets 3 mm × 7 mm FIRM, Cardiopoint™	88863213-56
	<b>2-0</b>	3	8×75	4 Blue and 4 White	6	Pledgets 3 mm × 3 mm FIRM, Cardiopoint™	88863203-56
	<b>2-0</b>	3	75	Blue	36	Cardiopoint™	88863202-51
	<b>2-0</b>	3	75	Blue	36	Cardiopoint™	88863088-51
	<b>2-0</b>	3	75	White	36	Cardiopoint™	88862855-51
	<b>2-0</b>	3	10×90	5 Blue and 5 White	6	Pledgets 3 mm × 7 mm FIRM, Cardiopoint™	88863290-56
	<b>2-0</b>	3	90	Blue	36	Pledgets 3 mm × 7 mm FIRM, Cardiopoint™	88863212-51
	<b>2-0</b>	3	90	Blue	36	Cardiopoint™	88863167-51
	<b>2-0</b>	3	90	White	36	Cardiopoint™	88862815-51
	<b>2-0</b>	3	10×90	5 Blue and 5 White	6	Cardiopoint™	88862852-56
	<b>2-0</b>	3	90	White	36	Pledgets 3 mm × 7 mm FIRM, Cardiopoint™	88862840-51
	<b>0</b>	3.5	8×75	4 Blue and 4 White	6	Pledgets 3 mm × 7 mm FIRM, Cardiopoint™	88863213-66
	<b>0</b>	3.5	75	Blue	36	Cardiopoint™	88863202-61
	<b>0</b>	3.5	90	Blue	36	Cardiopoint™	88863167-61





## TAPERCUTTING

Needles Specifications	Size		Length cm	Color	Units Per box	Comments	Reorder code	
	USP	EP						
<b>KV-15</b> 3/8 Circle 17 mm Double-arm		<b>2-0</b>	3	4×75	Blue	6	Surgalloy™	88863131-56
<b>KV-26</b> 3/8 Circle 26 mm Double-arm		<b>2-0</b>	3	10×90	5 Blue and 5 White	6	Pledgets 3 mm × 7 mm FIRM	88863275-56
<b>KV-5</b> 1/2 Circle 17 mm Double-arm		<b>4-0</b>	1.5	90	Blue	36	Surgalloy™	88863159-31
		<b>3-0</b>	2	90	Blue	36	Surgalloy™	88863159-41
		<b>2-0</b>	3	8×75	4 Blue and 4 White	6	Pledgets 3 mm × 3 mm FIRM Pledgets	88863236-56
		<b>2-0</b>	3	8×75	4 Blue and 4 White	6	Pledgets 3 mm × 7 mm FIRM	88863235-56
		<b>2-0</b>	3	10×75	5 Blue and 5 White	6	Surgalloy™	88862879-56
		<b>2-0</b>	3	90	Blue	36	Surgalloy™	88863159-51
		<b>2-0</b>	3	90	Blue	36	Pledgets 3 mm × 3 mm FIRM	88863193-51
<b>KV-7</b> 1/2 Circle 26 mm Double-arm		<b>2-0</b>	3	75	Blue	36	Surgalloy™	88863319-51
		<b>2-0</b>	3	8×75	4 Blue and 4 White	6	Pledgets 3 mm × 7 mm FIRM	88863240-56
		<b>2-0</b>	3	10×75	5 Blue and 5 White	6	Pledgets 3 mm × 7 mm FIRM	88863228-56
		<b>2-0</b>	3	90	Blue	36	Pledgets 3 mm × 7 mm FIRM	88863239-51
		<b>2-0</b>	3	90	Blue	36	Surgalloy™	88863168-51
<b>KV-34</b> 1/2 Circle 37 mm		<b>0</b>	3.5	75	Blue	36		88863119-61

## TAPERCUTTING

Needles Specifications	Size		Length cm	Color	Units Per box	Comments	Reorder code
	USP	EP					
<b>KV-37</b> 1/2 Circle    40 mm 	<b>0</b>	3.5	75	Blue	36		88863092-61
	<b>1</b>	4	75	Blue	36		88863092-71
	<b>2</b>	5	5×75	Blue	12		88863147-83
	<b>5</b>	7	5×75	Blue	12		88863147-89
<b>KV-40</b> 1/2 Circle    48 mm 	<b>5</b>	7	4×75	Blue	12		88862818-89

## STRAIGHT

Needles Specifications	Image	Size		Length cm	Color	Units Per box	Comments	Reorder code
		USP	EP					
<b>SS-24</b> 1/4 Circle 8.6 mm Double-arm Curvature 95°		6-0	0.7	45	White	12	Surgalloy™	CD-3113K
		5-0	1	45	White	12	Surgalloy™	CD-3123K
		5-0	1	45	White	12	Surgalloy™	D-1764K
		4-0	1.5	45	White	12	Surgalloy™	D-1763K
<b>SS-14</b> 1/4 Circle 8.6 mm Double-arm Curvature 115°		6-0	0.7	45	White	12	Surgalloy™	D-1722K
		5-0	1	45	White	12	Surgalloy™	D-1760K
<b>SS-28</b> 3/8 Circle 7.4 mm Double-arm Curvature 145°		6-0	0.7	45	White	12	Surgalloy™	D-1683K
<b>SS-2</b> 1/2 Circle 8.6 mm Double-arm Curvature 200°		5-0	1	45	White	12	Surgalloy™	CD-0423K
		5-0	1	45	White	12	Surgalloy™	D-1780K
		4-0	1.5	45	White	12	Surgalloy™	D-1779K

# VASCUFIL™

Coated Monofilament Polybutester

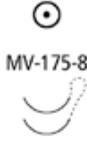


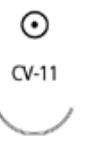



Innovation in Cardiovascular Surgery.






Please refer to the individual package IFUs for complete instructions, indications, contraindications, warnings and precautions.




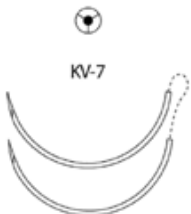


<b>Description</b>	Vascufil™ monofilament polybutester sutures are nonabsorbable, sterile, surgical sutures composed of a copolymer of butylene terephthalate and polytetramethylene ether glycol. The suture is coated with Polytribolate™, an absorbable polymer of e-caprolactone / glycolide / poloxamer 188.
<b>Indications</b>	Vascufil™ polybutester sutures are indicated for use in general soft tissue approximation and/or ligation, including use in cardiovascular and ophthalmic surgery, but not in microsurgery and neural tissue.
<b>Actions</b>	Vascufil™ sutures elicit a minimal acute inflammatory reaction in tissues, which is followed by gradual encapsulation of the suture by fibrous connective tissue.
<b>Tensile Strength</b>	Vascufil™ sutures are not absorbed, nor is any significant change in strength retention known to occur in vivo. Excretion of the absorbable coating is essentially complete within thirty days.
<b>Sterilization</b>	Gamma irradiation.
<b>Packaging</b>	Vascufil™ monofilament sutures are available dyed with Copper Phthalocyanine Blue in U.S.P. sizes 2–0 (3 Metric) through 7–0 (0.5 Metric). The sutures are supplied sterile, in pre-cut lengths affixed to various needle types, using permanent needle attachment techniques. Please refer to the individual package IFUs for complete instructions, indications, contraindications, warnings and precautions.

TAPER POINT								
Needles Specifications	Size USP	EP	Length cm	Color	Units Per box	Comments	Reorder code	
<b>MV-175-8</b> 3/8 Circle 8 mm Double-arm		5-0	1	75	Blue	36	Surgalloy™	88864741-21V
		7-0	0.5	4×60	Blue	12	Surgalloy™	88864788-05V
		7-0	0.7	2×60	Blue	12	Surgalloy™	88864795-06V
<b>CV-1</b> 3/8 Circle 9 mm Double-arm		7-0	0.5	45	Blue	36	Surgalloy™	88864703-01V
		7-0	0.5	60	Blue	36	Surgalloy™	88864705-01V
		7-0	0.5	4×60	Blue	12	Surgalloy™	88864702-05V
		7-0	0.5	75	Blue	36	Surgalloy™	88864704-01V
		7-0	0.5	4×75	Blue	12	Surgalloy™	88864701-05V
		6-0	0.7	45	Blue	36	Surgalloy™	88864703-11V
		6-0	0.7	60	Blue	36	Surgalloy™	88864705-11V
		6-0	0.7	75	Blue	36	Surgalloy™	88864704-11V
<b>CV</b> 3/8 Circle 11 mm Double-arm		6-0	0.7	75	Blue	36	Surgalloy™	88864756-11V
<b>CV-11</b> 3/8 Circle 13 mm		5-0	1	75	Blue	36	Surgalloy™	88864722-21V
<b>CV-11</b> 3/8 Circle 13 mm Double-arm		6-0	0.7	45	Blue	36	Surgalloy™	88864707-11V
		6-0	0.7	60	Blue	36	Surgalloy™	88864709-11V
		6-0	0.7	75	Blue	36	Surgalloy™	88864710-11V
		5-0	1	60	Blue	36	Surgalloy™	88864709-21V
		5-0	1	75	Blue	36	Surgalloy™	88864710-21V
		5-0	1	90	Blue	36	Surgalloy™	88864711-21V
<b>CV-15</b> 3/8 Circle 17 mm		5-0	1	75	Blue	36	Surgalloy™	88864739-21V
		4-0	1.5	75	Blue	36	Surgalloy™	88864739-31V
<b>CV-15</b> 3/8 Circle 17 mm Double-arm		5-0	1	90	Blue	36	Surgalloy™	88864716-21V
		4-0	1.5	90	Blue	36	Surgalloy™	88864716-31V

## TAPER POINT

Needles Specifications	Size USP	EP	Length cm	Color	Units Per box	Comments	Reorder code	
<b>CV-22</b> 1/2 Circle 13 mm Double-arm		6-0	0.7	4×75	Blue	12	Surgalloy™	88864742-15V
		6-0	0.7	75	Blue	36	Surgalloy™	88864741-11V
		7-0	0.5	4×60	Blue	36	Surgalloy™	88864794-01V
<b>CV-23</b> 1/2 Circle 17 mm Double-arm		5-0	1	90	Blue	36	Surgalloy™	88864719-21V
		4-0	1.5	90	Blue	36	Surgalloy™	88864719-31V
<b>V-20</b> 1/2 Circle 26 mm Double-arm		2-0	3	90	Blue	36	Surgalloy™	88864750-51V
		3-0	2	90	Blue	36	Surgalloy™	88864750-41V
		4-0	1.5	90	Blue	36	Surgalloy™	88864750-31V

## TAPERCUTTING

Needles Specifications	Size USP	EP	Length cm	Color	Units Per box	Comments	Reorder code	
<b>KV-7</b> 1/2 Circle 26 mm		2-0	3	90	Blue	36	Surgalloy™	88864730-51V
<b>KV-7</b> 1/2 Circle 26 mm Double-arm		4-0	1.5	90	Blue	36	Surgalloy™	88864731-31V
		3-0	2	90	Blue	36	Surgalloy™	88864731-41V
		2-0	3	90	Blue	36	Surgalloy™	88864731-51V
<b>KV-11</b> 3/8 Circle 13 mm Double-arm		6-0	0.7	75	Blue	36	Surgalloy™	88864734-11V
		5-0	1	75	Blue	36	Surgalloy™	88864734-21V
<b>KV-1</b> 3/8 Circle 9 mm Double-arm		6-0	0.7	60	Blue	36	Surgalloy™	88864714-11V

# FLEXON™

## Multifilament Temporary Cardiac Pacing Lead

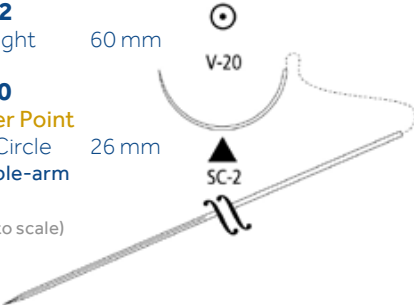
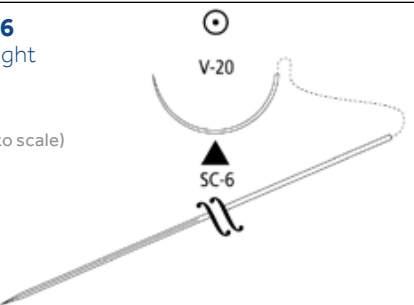
Excellent Stability with Great Flexibility.



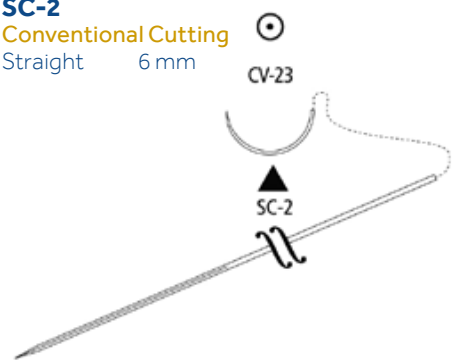
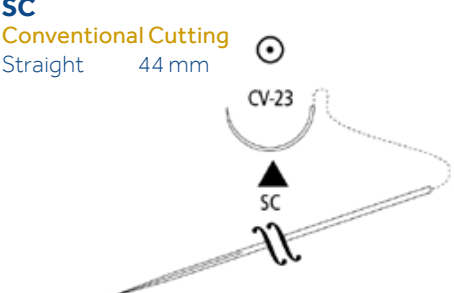
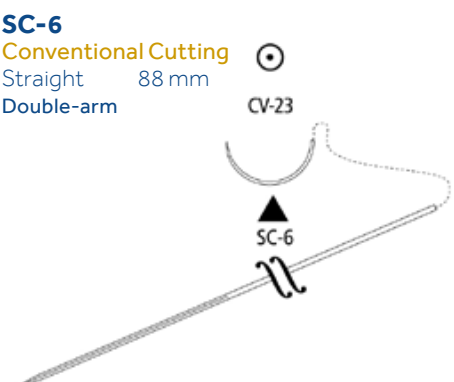
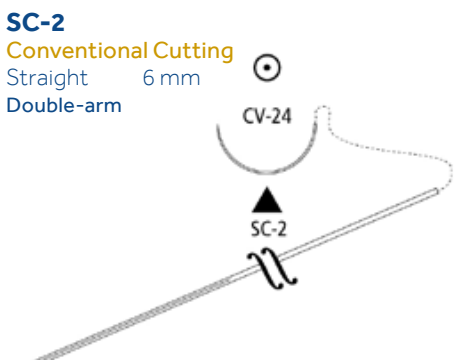
Please refer to the individual package IFUs for complete instructions, indications, contraindications, warnings and precautions.

<b>Description</b>	Flexon™ temporary cardiac pacing lead is a stainless steel twisted, multistrand wire coated with orange or white PTFE poly (tetrafluoroproethylene) or clear FEP poly (tetrafluoroethylene-co- hexafluoropropylene).
<b>Indications</b>	Flexon™ temporary cardiac pacing lead is indicated for use in temporary atrial and ventricular pacing and sensing during and after cardiac surgery.
<b>Actions</b>	Flexon™ temporary cardiac pacing lead elicits minimal acute inflammatory reaction in tissues and is not absorbed.
<b>Tensile Strength</b>	Permanent.
<b>Sterilization</b>	Ethylene oxide.
<b>Packaging</b>	Flexon™ temporary cardiac pacing lead are available coated with orange or white PTFE or clear FEP in size 0 (3.5 metric), size 2-0 (3 metric) and 3-0 (2 metric). Flexon™ temporary cardiac pacing lead is supplied sterile affixed to various types of needles, including Breakaway™ needles.

TAPER POINT

Needles Specifications	Size		Length cm	Color	Units Per box	Comments	Reorder code
	USP	EP					
<p><b>SC-2</b> Straight 60 mm</p> <p><b>V-20</b> Taper Point 1/2 Circle 26 mm Double-arm</p> <p>(Not to scale)</p> 	<b>0</b>	3.5	60	Undyed	12		88862597-63
<p><b>SC-6</b> Straight</p> <p>(Not to scale)</p> 	<b>2-0</b>	3	2 × 60	1 Orange 1 White	12		88862623-53
	<b>2-0</b>	3	60	Orange	12		88862617-53
	<b>2-0</b>	3	60	White	12		88862615-53
	<b>0</b>	3.5	60	Undyed	12		88862589-63

## TAPER POINT

Needles Specifications	Size		Length cm	Color	Units Per box	Comments	Reorder code
	USP	EP					
<b>CV-23</b> 1/2 Circle 17 mm  <b>SC-2</b> Conventional Cutting Straight 6 mm 	3-0	2	60	Undyed	12	FEP Coating	88862592-43
<b>SC</b> Conventional Cutting Straight 44 mm 		2	60	Undyed	12	FEP Coating	88862591-43
<b>SC-6</b> Conventional Cutting Straight 88 mm Double-arm 		2	60	Undyed	12	Breakaway™ FEP Coating	88862586-43
<b>CV-24</b> 1/2 Circle 20 mm  <b>SC-2</b> Conventional Cutting Straight 6 mm Double-arm 	0	3.5	60	Undyed	12	FEP Coating	88862593-63

# STEEL

## Monofilament Stainless Steel

Providing Maximum Tensile Strength.



Please refer to the individual package IFUs for complete instructions, indications, contraindications, warnings and precautions.



**Description**

Stainless monofilament steel sutures are composed of 316L stainless steel conforming to ASTM Standard F138 Grade 2 "Stainless steel bar and wire for surgical implants."

**Indications**

Steel sutures are indicated for use in abdominal wound closure, hernia repair and sternal closure.

**Actions**

Steel sutures elicit a minimal acute inflammatory reaction in tissue, and are not absorbed.

**Tensile Strength**

Permanent.

**Sterilization**

Ethylene oxide or gamma irradiation.

**Packaging**


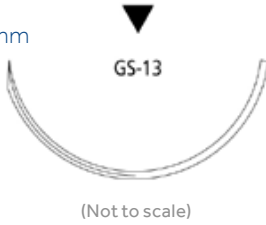
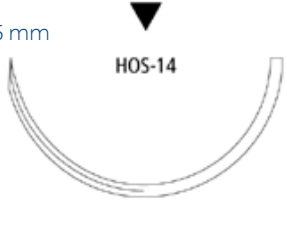
Steel sutures are available in monofilament strands in size 1 (4 Metric), size 2 (6 metric), size 4 (6 metric) through 7 (9 Metric). The sutures are supplied sterile, in pre-cut lengths affixed to various types of needles using permanent needle attachment technique.

# Steel

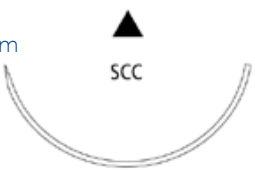
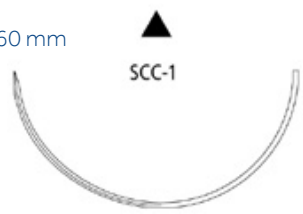
## PRE-CUT

Needles Specifications	Size		Length cm	Color	Units Per box	Comments	Reorder code
	USP	EP					
	<b>2-0</b>	3	12 × 45	Undyed	12	Pre-cut	88862407-53
	<b>0</b>	3.5	12 × 45	Undyed	12	Pre-cut	88862407-63
	<b>1</b>	4	6 × 45	Undyed	12	Pre-cut	88862408-73
	<b>2</b>	5	6 × 45	Undyed	12	Pre-cut	88862408-83
	<b>3-0</b>	2	12 × 45	Undyed	12	Pre-cut	88862407-43
	<b>7</b>	9	12 × 45	Undyed	12	Pre-cut	88862408-73


## REVERSE CUTTING

Needles Specifications	Size		Length cm	Color	Units Per box	Comments	Reorder code
	USP	EP					
<b>GS-12</b> 1/2 Circle 40 mm 	<b>2</b>	5	5 × 45	Undyed	12		88862413-83
<b>GS-13 ROTO</b> 1/2 Circle 48 mm  (Not to scale)	<b>2</b>	5	2 × 45	Undyed	12	ROTOGRIP™	88862420-83
	<b>4</b>	6	2 × 45	Undyed	12	ROTOGRIP™	88862420-09
<b>HOS-14</b> 1/2 Circle 55 mm  (Not to scale)	<b>5</b>	7	4 × 45	Undyed	12	ROTOGRIP™	88862410-89
	<b>5</b>	7	2 × 75	Undyed	12	ROTOGRIP™	88862424-89
	<b>5</b>	7	75	Undyed	12	ROTOGRIP™	88862412-89
	<b>6</b>	8	4 × 45	Undyed	12	ROTOGRIP™	88862410-69



## CONVENTIONAL CUTTING

Needles Specifications	Size		Length cm	Color	Units Per box	Comments	Reorder code
	USP	EP					
<b>SCC</b> 1/2 Circle    48 mm 	5	7	4 × 45	Undyed	12		88862224-89
	5	7	2 × 45	Undyed	12		88862222-89
	5	7	6 × 45	Undyed	12		88862226-89
	6	8	4 × 45	Undyed	12		88862224-69
	6	8	6 × 45	Undyed	12		88862226-69
	6	8	2 × 45	Undyed	12		88862222-69
	7	9	4 × 45	Undyed	12		88862224-49
<b>SCC-1</b> 1/2 Circle    60 mm 	5	7	2 × 75	Undyed	12		88862228-89

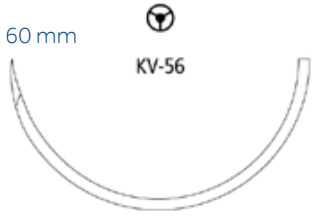
## TAPER POINT

Needles Specifications	Size		Length cm	Color	Units Per box	Comments	Reorder code
	USP	EP					
<b>GS-21</b> 1/2 Circle 37 mm 	2-0	3	5 × 45	Undyed	12		88862345-53
<b>GS-25</b> 1/2 Circle 48 mm 	2-0	3	5 × 60	Undyed	12		88862348-53



## TAPERCUTTING

Needles Specifications	Size		Length cm	Color	Units Per box	Comments	Reorder code
	USP	EP					
<b>KV-37</b> 1/2 Circle 40 mm 	2	5	5 × 45	Undyed	12		88862382-83
<b>KV-40</b> 1/2 Circle 48 mm  (Not to scale)	3	6	5 × 45	Undyed	12		88862466-25
	5	7	4 × 45	Undyed	12		88862414-89
	5	7	4 × 45	Undyed	12	ROTOGRIP	88862393-89
	5	7	2 × 45	Undyed	12		88862376-89
	6	8	2 × 45	Undyed	12		88862392-69
	6	8	4 × 45	Undyed	12		88862393-69
	7	9	2 × 45	Undyed	12		88862392-49


## TAPERCUTTING

Needles Specifications	Size		Length cm	Color	Units Per box	Comments	Reorder code
	USP	EP					
<b>KV-56</b> 1/2 Circle    60 mm 	5	7	5 × 75	Undyed	12		88862396-89

## BLUNT TAPER POINT

Needles Specifications	Size		Length cm	Color	Units Per box	Comments	Reorder code
	USP	EP					
<b>BGS-26</b> 1/2 Circle    44 mm  (Not to scale)	5	7	75	Undyed	12	ProtectPoint™	88862497-89
<b>BGS-29</b> 1/2 Circle    63 mm 	5	7	5 × 45	Undyed	12	ProtectPoint™	88862495-89

## BLUNT POINT

Needles Specifications	Size		Length cm	Color	Units Per box	Comments	Reorder code
	USP	EP					
<b>BP-27</b> 1/2 Circle    64 mm 	5	7	2 × 75	None	12	ProtectPoint™	88862494-89

# DERMALON™ / MONOSOF™

Monofilament Nylon

High Tensile Strength. Low tissue Reactivity.




Please refer to the individual package IFUs for complete instructions, indications, contraindications, warnings and precautions.

<b>Description</b>	DERMALON™ and MONOSOF™ monofilament sutures are inert, nonabsorbable, sterile surgical sutures composed of the longchain aliphatic polymers Nylon 6 and Nylon 6.6.
<b>Indications</b>	DERMALON™ and MONOSOF™ monofilament nylon sutures are indicated for use in general soft tissue approximation and/or ligation, including use in cardiovascular, ophthalmic, and neurological surgery.
<b>Actions</b>	Both DERMALON™ and MONOSOF™ monofilament nylon sutures elicit a minimal acute inflammatory reaction in tissue, which is followed by gradual encapsulation of the suture by fibrous connective tissue.
<b>Tensile Strength</b>	As with any nylon suture, though not absorbed, progressive hydrolysis of the suture, in vivo, may result in gradual loss of its tensile strength over time.
<b>Sterilization</b>	DERMALON™ and MONOSOF™ monofilament nylon sutures are sterilized by ethylene oxide or gamma irradiation.
<b>Packaging</b>	DERMALON™ monofilament nylon sutures are available dyed blue with FD&C Blue No. 2 in sizes 2-0 (3 Metric) through 6-0 (0.7 metric). MONOSOF™ monofilament nylon sutures are available either clear (undyed) or dyed black with Logwood extract, in sizes 2 (5 Metric) through 11-0 (0.1 Metric). Both sutures are supplied sterile, in pre-cut lengths affixed to needles using permanent needle attachment techniques.

## Dermalon™ / Monosof™



Dermalon™ / Monosof™ suture is a brand you have trusted for ~ 50 years with a proven history of performance.

PRE-CUT							
Needles Specifications	Size		Length cm	Color	Units Per box	Comments	Reorder code
	USP	EP					
	<b>3-0</b>	2	12×45	Black	36	Pre-cut	N-63
	<b>2-0</b>	3	6×45	Black	36	Pre-cut	N-59
	<b>2</b>	5	6×45	Black	36	Pre-cut	N-64




REVERSE CUTTING							
Needles Specifications	Size		Length cm	Color	Units Per box	Comments	Reorder code
	USP	EP					
<b>C-1</b> 3/8 Circle 12 mm 	<b>6-0</b>	0.7	25	Black	12		SN-770G
	<b>6-0</b>	0.7	45	Black	12		SN-667G
	<b>6-0</b>	0.7	45	Black	36		SN-667
	<b>6-0</b>	0.7	45	Blue	36		88861741-11
	<b>6-0</b>	0.7	75	Black	36		SN-367
	<b>5-0</b>	1	25	Black	12		SN-771G
	<b>5-0</b>	1	45	Black	36		SN-871
	<b>5-0</b>	1	45	Blue	36		88861741-21
	<b>5-0</b>	1	75	Black	36		SN-368
	<b>4-0</b>	1.5	25	Black	12		SN-772G
	<b>4-0</b>	1.5	45	Black	36		SN-872
	<b>4-0</b>	1.5	45	Blue	36		88861741-31




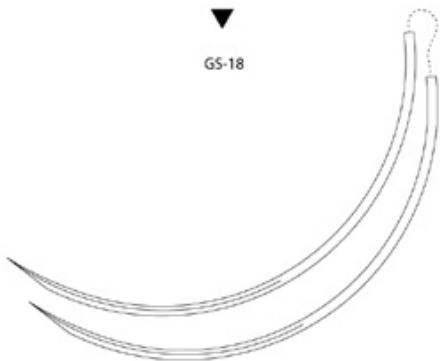
## REVERSE CUTTING

Needles Specifications	Size		Length cm	Color	Units Per box	Comments	Reorder code
	USP	EP					
<b>C-12</b> 3/8 Circle 16 mm 	6-0	0.7	45	Blue	36		88861747-11
	6-0	0.7	45	Black	36		SN-660
	6-0	0.7	75	Black	36		SN-760
	5-0	1	45	Black	36		SN-659
	5-0	1	45	Black	12		SN-659G
	5-0	1	75	Black	36		SN-761
	4-0	1.5	75	Black	36		SN-762
	3-0	2	75	Black	36		SN-763
<b>C-13</b> 3/8 Circle 19 mm 	6-0	0.7	45	Blue	36		88861756-11
	5-0	1	45	Blue	36		88861756-21
	5-0	1	45	Black	12		SN-661G
	5-0	1	45	Black	36		SN-661
	5-0	1	75	Black	36		SN-643
	4-0	1.5	25	Black	12		SN-552G
	4-0	1.5	45	Black	12		SN-662G
	4-0	1.5	45	Black	36		SN-662
	4-0	1.5	45	Blue	36		88861756-31
	4-0	1.5	75	Blue	36		88861757-31
	4-0	1.5	90	Black	36		SN-644
	3-0	2	45	Black	36		SN-653
	3-0	2	45	Blue	36		88861756-41
	3-0	2	75	Black	36		SN-641
	3-0	2	75	Blue	36		88861757-41
	3-0	2	90	Black	36		SN-665
	2-0	3	45	Blue	36		88861756-51
	2-0	3	75	Black	36		SN-750
	2-0	3	90	Black	36		SN-666




## REVERSE CUTTING

Needles Specifications	Size		Length cm	Color	Units Per box	Comments	Reorder code
	USP	EP					
<b>C-14</b> 3/8 Circle 24 mm 	4-0	1.5	45	Black	36		SN-658
	4-0	1.5	45	Blue	36		88861744-31
	4-0	1.5	75	Blue	36		88861799-31
	4-0	1.5	90	Black	36		SN-649
	3-0	2	45	Black	12		SN-663G
	3-0	2	45	Black	36		SN-663
	3-0	2	45	Blue	36		88861744-41
	3-0	2	75	Blue	36		88861799-41
	3-0	2	75	Black	36		SN-669
	3-0	2	90	Black	36		SN-648
	2-0	3	45	Black	36		SN-758
	2-0	3	75	Blue	36		88861799-51
	2-0	3	90	Black	36		SN-747
<b>C-15</b> 3/8 Circle 26 mm 	2-0	3	45	Black	12		SN-764G
	2-0	3	45	Black	36		SN-764
	2-0	3	45	Blue	36		88861744-51
	2-0	3	75	Black	12		SN-664G
	2-0	3	75	Black	36		SN-664
	2-0	3	90	Blue	36		88861720-51
	0	3.5	75	Black	36		SN-765
<b>C-16</b> 3/8 Circle 30 mm 	4-0	1.5	75	Black	36		SN-670
	3-0	2	75	Blue	36		88861750-41
	3-0	2	75	Black	12		SN-671G
	3-0	2	75	Black	36		SN-671
	3-0	2	90	Black	36		SN-771
	2-0	3	75	Black	36		SN-688
	2-0	3	90	Black	36		SN-788
	0	3.5	75	Black	36		SN-689
	1	4	90	Black	36		SN-790





## REVERSE CUTTING

Needles Specifications	Size		Length cm	Color	Units Per box	Comments	Reorder code
	USP	EP					
<b>C-18</b> 3/8 Circle 35 mm 	<b>3-0</b>	2	90	Black	36		SN-693
	<b>2-0</b>	3	90	Black	36		SN-694
<b>C-17</b> 3/8 Circle 39 mm 	<b>4-0</b>	1.5	75	Black	36		SN-672
	<b>3-0</b>	2	75	Black	12		SN-673G
	<b>3-0</b>	2	75	Black	36		SN-673
	<b>3-0</b>	2	75	Blue	36		88861727-41
	<b>2-0</b>	3	75	Blue	36		88861727-51
	<b>2-0</b>	3	75	Black	36		SN-674
	<b>0</b>	3.5	90	Black	36		SN-775
	<b>1</b>	4	90	Black	36		SN-776
	<b>2</b>	5	90	Black	36		SN-777
<b>GS-18</b> 3/8 Circle 77 mm 	<b>2</b>	5	45	Black	24	With bolster and weight	CN-290
<b>GS-18</b> 3/8 Circle 77 mm Double-arm 	<b>2</b>	5	75	Black	24		CN-490


## REVERSE CUTTING

Needles Specifications	Image	Size		Length cm	Color	Units Per box	Comments	Reorder code
		USP	EP					
<b>C-22</b> 1/2 Circle    18 mm		5-0	1	75	Black	36		SN-629
		4-0	1.5	75	Black	36		SN-630
		3-0	2	75	Black	36		SN-631
<b>C-23</b> 1/2 Circle    24 mm		3-0	2	90	Black	36		SN-247
<b>C-25</b> 1/2 Circle    30 mm		3-0	2	75	Black	36		SN-681
		3-0	2	90	Black	36		SN-781
		2-0	3	75	Black	36		SN-682
		2-0	3	90	Black	36		SN-782
<b>GS-11</b> 1/2 Circle    37 mm		3-0	2	90	Black	36		SN-881
		2-0	3	75	Black	36		CN-819
		1	4	90	Black	36		CN-822



## PREMIUM REVERSE CUTTING

Needles Specifications	Size USP	EP	Length cm	Color	Units Per box	Comments	Reorder code
<b>P-16</b> 3/8 Circle 7 mm 	7-0	0.5	45	Black	12		SN-1647G
	7-0	0.5	45	Black	36		SN-1647
<b>P-10</b> 3/8 Circle 11 mm 	7-0	0.5	45	Black	12	Surgalloy™	SN-3686
	6-0	0.7	45	Black	12	Surgalloy™	SN-3697
	6-0	0.7	45	Clear	12	Surgalloy™	SN-1689G
	6-0	0.7	45	Black	36	Surgalloy™	SN-5697
	6-0	0.7	45	Clear	36	Surgalloy™	SN-1689
	6-0	0.7	45	Blue	12		88861754-13
	5-0	1	45	Black	12	Surgalloy™	SN-3695
<b>P-13</b> 3/8 Circle 13 mm 	6-0	0.7	45	Black	36	Surgalloy™	SN-5696
	6-0	0.7	45	Black	12	Surgalloy™	SN-5696G
	6-0	0.7	45	Blue	12	Surgalloy™	SDN-5696G
	5-0	1	45	Black	12	Surgalloy™	SN-5698G
	5-0	1	45	Clear	12	Surgalloy™	SN-5690G
	5-0	1	45	Black	36	Surgalloy™	SN-5698
	5-0	1	45	Clear	36	Surgalloy™	SN-5690
	5-0	1	45	Blue	12	Surgalloy™	88861749-24
	4-0	1.5	45	Black	12	Surgalloy™	SN-5699G
	4-0	1.5	45	Black	36	Surgalloy™	SN-5699
	4-0	1.5	45	Clear	36	Surgalloy™	SN-5691
	4-0	1.5	45	Blue	12	Surgalloy™	SDN-5691G
<b>P-11</b> 3/8 Circle 16 mm 	6-0	0.7	45	Black	12	Surgalloy™	SN-5767G
	5-0	1	45	Black	12	Surgalloy™	SN-5668G
	5-0	1	45	Black	36	Surgalloy™	SN-5668
	4-0	1.5	45	Black	12	Surgalloy™	SN-5694G
	4-0	1.5	45	Black	36	Surgalloy™	SN-5694

## PREMIUM REVERSE CUTTING

Needles Specifications	Size		Length cm	Color	Units Per box	Comments	Reorder code
	USP	EP					
<b>P-12</b> 3/8 Circle    19 mm 	6-0	0.7	45	Black	12	Surgalloy™	SN-5665G
	5-0	1	45	Black	12	Surgalloy™	SN-5666G
	5-0	1	45	Black	36	Surgalloy™	SN-5666
	4-0	1.5	45	Clear	12	Surgalloy™	SN-5677
	4-0	1.5	45	Black	12	Surgalloy™	SN-5667G
	4-0	1.5	45	Black	36	Surgalloy™	SN-5667
	4-0	1.5	45	Blue	12	Surgalloy™	SDN-5666G
	3-0	2	45	Black	12	Surgalloy™	SN-5669G
	3-0	2	45	Black	36	Surgalloy™	SN-5669
	2-0	3	45	Black	36	Surgalloy™	SN-5670

## PREMIUM REVERSE CUTTING



Needles Specifications	Size		Length cm	Color	Units Per box	Comments	Reorder code
	USP	EP					
<b>P-14</b> 3/8 Circle    24 mm 	<b>4-0</b>	1.5	45	Black	12	Surgalloy™	SN-5862G
	<b>3-0</b>	2	45	Black	12	Surgalloy™	SN-5663G
	<b>3-0</b>	2	45	Black	36	Surgalloy™	SN-5663
	<b>3-0</b>	2	45	Blue	12	Surgalloy™	SDN-5943G
<b>P-24</b> 1/2 Circle    16 mm 	<b>6-0</b>	0.7	45	Black	12	Surgalloy™	SN-5660G
	<b>4-0</b>	1.5	45	Black	12	Surgalloy™	SN-5662G
	<b>4-0</b>	1.5	45	Black	36	Surgalloy™	SN-5662

## CONVENTIONAL CUTTING















Needles Specifications	Size USP	EP	Length cm	Color	Units Per box	Comments	Reorder code
<b>PC-13</b> 3/8 Circle 11 mm 	<b>6-0</b>	0.7	45	Black	36		SN-1956
	<b>5-0</b>	1	45	Black	36		SN-1955
	<b>4-0</b>	1.5	45	Black	36		SN-1954
<b>PC-11</b> 3/8 Circle 16 mm 	<b>6-0</b>	0.7	45	Black	36		SN-1966
	<b>5-0</b>	1	45	Black	12	Surgalloy™	SN-3965
	<b>4-0</b>	1.5	45	Black	36		SN-1964
<b>PC-12</b> 3/8 Circle 19 mm 	<b>5-0</b>	1	45	Black	36		SN-1995
	<b>4-0</b>	1.5	45	Black	36		SN-1994
	<b>3-0</b>	2	45	Black	36		SN-1693
<b>SCE-4</b> 3/8 Circle 19 mm 	<b>5-0</b>	1	45	Blue	36		88861742-21
	<b>4-0</b>	1.5	45	Blue	36		88861742-31
	<b>4-0</b>	1.5	45	Black	36		SN-568
	<b>3-0</b>	2	45	Blue	36		88861742-41
<b>GCC-90</b> 3/8 Circle 90 mm 	<b>2-0</b>	3	100	Black	36		CN-791
	<b>0</b>	3.5	100	Black	36		CN-792
	<b>1</b>	4	100	Black	36		CN-793



## CONVENTIONAL CUTTING

Needles Specifications		Size		Length cm	Color	Units Per box	Comments	Reorder code
		USP	EP					
<b>SC-1</b> Straight      51 mm 		<b>4-0</b>	1.5	90	Black	36		SN-615
		<b>3-0</b>	2	90	Black	36		SN-616
		<b>2-0</b>	3	90	Black	36		SN-617
<b>SC-2</b> Straight      60 mm 		<b>3-0</b>	2	75	Black	36		SN-627
		<b>3-0</b>	2	75	Blue	36		88861764-41
		<b>2-0</b>	3	75	Black	36		SN-628
		<b>2-0</b>	3	75	Blue	36		88861764-51
		<b>2-0</b>	3	100	Black	36		SN-728


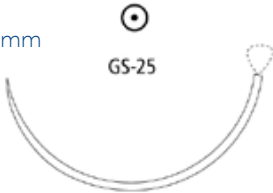
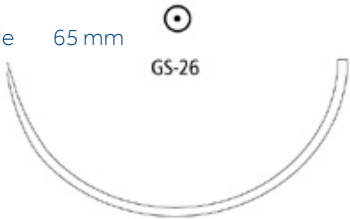

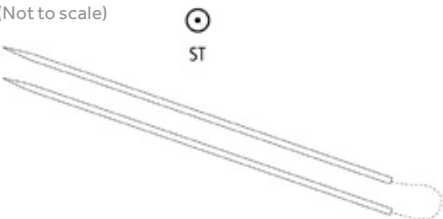
## TAPER POINT

Needles Specifications			Size		Length cm	Color	Units Per box	Comments	Reorder code
			USP	EP					
<b>MV-50-3</b> 3/8 Circle 4 mm	 MV-50-3 	<b>11-0</b>	0.1	13	Black	12	Surgalloy™	N-2500	
<b>MV-70-3</b> 3/8 Circle 4 mm	 MV-70-3 	<b>11-0</b>	0.1	13	Black	12	Surgalloy™	N-2510	
		<b>10-0</b>	0.2	13	Black	12	Surgalloy™	N-2512	
<b>MV-135-3</b> 3/8 Circle 4 mm	 MV-135-3 	<b>10-0</b>	0.2	13	Black	12	Surgalloy™	N-2540	
		<b>9-0</b>	0.3	8	Black	12	Surgalloy™	N-1603G	
		<b>8-0</b>	0.4	13	Black	12	Surgalloy™	N-2541	
<b>MV-70-4</b> 3/8 Circle 5 mm	 MV-70-4 	<b>10-0</b>	0.2	13	Black	12	Surgalloy™	N-2514	
<b>MV-100-4</b> 3/8 Circle 5 mm	 MV-100-4 	<b>10-0</b>	0.2	13	Black	12	Surgalloy™	N-2530	
		<b>9-0</b>	0.3	13	Black	12	Surgalloy™	N-2532	
<b>MV-135-4</b> 3/8 Circle 5 mm	 MV-135-4 	<b>10-0</b>	0.2	13	Black	12	Surgalloy™	N-2544	
		<b>9-0</b>	0.3	13	Black	12	Surgalloy™	N-2546	
		<b>8-0</b>	0.4	13	Black	12	Surgalloy™	N-2543	
<b>MV-135-5</b> 3/8 Circle 6 mm	 MV-135-5 	<b>9-0</b>	0.3	13	Black	12	Surgalloy™	N-2548	
		<b>8-0</b>	0.4	13	Black	12	Surgalloy™	N-2547	


## TAPER POINT

Needles Specifications	Size		Length cm	Color	Units Per box	Comments	Reorder code
	USP	EP					
<b>GS-22</b> 1/2 Circle 27 mm 	0	3.5	100	Black	36		CN-225
<b>GS-24</b> 1/2 Circle 40 mm 	0	3.5	100	Black	36		CN-723
	1	4	100	Black	36		CN-722
<b>GS-24 LOOP</b> 1/2 Circle 40 mm 	0	3.5	150	Black	24	Loop	CN-825-L



## TAPER POINT

Needles Specifications	Size		Length cm	Color	Units Per box	Comments	Reorder code
	USP	EP					
<b>GS-25</b> 1/2 Circle 48 mm 	1	4	100	Black	36		CN-724
<b>GS-25 LOOP</b> 1/2 Circle 48 mm 	1	4	150	Black	24	Loop	CN-830-L
<b>GS-26</b> 1/2 Circle 65 mm 	1	4	150	Black	24		CN-824
	2	5	150	Black	24		CN-825
<b>GS-26 LOOP</b> 1/2 Circle 65 mm 	0	3.5	240	Black	24	Loop	CN-823-L
<b>ST</b> Straight 60 mm Double-arm (Not to scale) 	2-0	3	90	Black	36		GN-283
	0	3.5	90	Black	36		GN-284


## HEAVY TAPER POINT

Needles Specifications	Size		Length cm	Color	Units Per box	Comments	Reorder code
	USP	EP					
<b>HGS-22</b> 1/2 Circle    27 mm 	0	3.5	100	Black	36		CN-925
	1	4	100	Black	36		CN-926

## TAPERCUTTING

Needles Specifications	Size		Length cm	Color	Units Per box	Comments	Reorder code
	USP	EP					
<b>MVK-70-3</b> 3/8 Circle    4 mm 	10-0	0.2	13	Black	12	Surgalloy™	N-2513
<b>MVK-100-4</b> 3/8 Circle    5 mm 	10-0	0.2	13	Black	12	Surgalloy™	N-2531
	9-0	0.3	13	Black	12	Surgalloy™	N-2533

## BLUNT POINT

Needles Specifications	Size		Length cm	Color	Units Per box	Comments	Reorder code
	USP	EP					
<b>BGS-25</b> 1/2 Circle    48 mm 	1	4	100	Black	36	ProtectPoint™	CN-771

## PREMIUM SPATULA

Needles Specifications	Image	Size		Length cm	Color	Units Per box	Comments	Reorder code
		USP	EP					
<b>SE-100-8</b> 1/4 Circle 6.3 mm Double-arm Curvature 100°		8-0	0.4	30	Black	12	Surgalloy™	N-2714K
<b>SS-24</b> 1/4 Circle 8.6 mm Double-arm Curvature 95°		5-0	1	45	Black	12	Surgalloy™	N-2781K
<b>SE-140-6</b> 3/8 Circle 6.3 mm Curvature 140°		10-0	0.2	20	Black	12	Surgalloy™	N-2730-K
<b>SE-140-6</b> 3/8 Circle 6.3 mm Double-arm Curvature 140°		10-0	0.2	30	Black	12	Surgalloy™	N-2770-K
		9-0	0.3	30	Black	12	Surgalloy™	N-2771-K
<b>SE-140-8</b> 3/8 Circle 6.3 mm Curvature 140°		10-0	0.2	20	Black	12	Surgalloy™	N-2736K
<b>SE-140-8</b> 3/8 Circle 6.3 mm Double-arm Curvature 140°		10-0	0.2	30	Black	12	Surgalloy™	N-2718K
		9-0	0.3	30	Black	12	Surgalloy™	N-2717-K
<b>SE-CC-6</b> Compound curve 4.8 mm Double-arm Curvature 130°		10-0	0.2	30	Black	12	Surgalloy™	N-2799-K
<b>SE-160-4</b> 1/2 Circle 5.5 mm Double-arm Curvature 160°		11-0	0.1	30	Black	12	Surgalloy™	N-2705K
		10-0	0.2	30	Black	12	Surgalloy™	N-2704K

## SPATULA

Needles Specifications	Image	Size		Length cm	Color	Units Per box	Comments	Reorder code
		USP	EP					
<b>SE-160-6</b> 1/2 Circle 5.5 mm Curvature 160°	 SE-160-6 	10-0	0.2	20	Black	12	Surgalloy™	N-2734K
<b>SE-160-6</b> 1/2 Circle 5.5 mm Double-arm Curvature 160°	 SE-160-6  	10-0	0.2	15	Black	12	Surgalloy™	N-2757K
		10-0	0.2	30	Black	12	Surgalloy™	N-2756K
		9-0	0.3	30	Black	12	Surgalloy™	N-2760K
<b>SE-160-8</b> 1/2 Circle 5.5 mm Curvature 160°	 SE-160-8 	10-0	0.2	15	Black	12	Surgalloy™	N-2766K
<b>SE-160-8</b> 1/2 Circle 5.5 mm Double-arm Curvature 160°	 SE-160-8  	10-0	0.2	30	Black	12	Surgalloy™	N-2786K
		9-0	0.3	30	Black	12	Surgalloy™	N-2787K
<b>SE-175-8</b> 1/2 Circle 6.9 mm Double-arm Curvature 175°	 SE-175-8  	10-0	0.2	30	Black	12	Surgalloy™	N-2719K
		8-0	0.4	30	Black	12	Surgalloy™	N-2716K
<b>SE-175-6</b> 1/2 Circle 7 mm Double-arm Curvature 175°	 SE-175-6  	10-0	0.2	30	Black	12		N-2966K
		9-0	0.3	30	Black	12	Surgalloy™	N-2967
<b>SS-2</b> 1/2 Circle 8.6 mm Double-arm Curvature 200°	 SS-2  	5-0	1	45	Black	12	Surgalloy™	N-2761K

# NOVAFIL™

Monofilament Polybutester

Supple and Flexible. Stretches with Healing.







Please refer to the individual package IFUs for complete instructions, indications, contraindications, warnings and precautions.







<b>Description</b>	Novafil™ monofilament polybutester sutures are nonabsorbable, sterile, surgical sutures composed of a copolymer of butylene terephthalate and polytetramethylene ether glycol.
<b>Indications</b>	Novafil™ polybutester sutures are indicated for use in general soft tissue approximation and/or ligation, including use in cardiovascular and ophthalmic surgery, but not in microsurgery and neural tissue.
<b>Actions</b>	Novafil™ polybutester sutures elicit a minimal acute inflammatory reaction in tissues, followed by gradual encapsulation of the suture by fibrous connective tissue.
<b>Tensile Strength</b>	Novafil™ polybutester sutures are not absorbed, nor is any significant change in tensile strength retention known to occur in vivo.
<b>Sterilization</b>	Gamma irradiation.
<b>Packaging</b>	Novafil™ polybutester sutures are available clear or pigmented with Copper Phthalocyanine Blue, in sizes 1 (4 Metric), size 2-0 (metric 3) through 7-0 (0.5 Metric). The sutures are supplied sterile, in pre-cut lengths affixed to reverse cutting needle type using permanent needle attachment techniques. Please refer to the individual package IFUs for complete instructions, indications, contraindications, warnings and precautions.

# Novafil™





Novafil™ suture is a brand you have trusted for ~ 30 years with a proven history of performance.

REVERSE CUTTING								
Needles Specifications	Size USP	EP	Length cm	Color	Units Per box	Comments	Reorder code	
<b>C-1</b> 3/8 Circle 12 mm 	<b>6-0</b>	0.7	45	Blue	12		8886441913	
	<b>5-0</b>	1	45	Blue	12		8886441923	
	<b>4-0</b>	1.5	45	Blue	12		8886441933	
<b>C-12</b> 3/8 Circle 16 mm 	<b>5-0</b>	1	45	Blue	12		88864555-23	
<b>C-13</b> 3/8 Circle 19 mm 	<b>6-0</b>	0.7	45	Blue	12	455	8886442013	
	<b>5-0</b>	1	45	Blue	12		8886442023	
	<b>4-0</b>	1.5	45	Blue	12		8886442033	
	<b>4-0</b>	1.5	75	Blue	36		8886403231	
	<b>3-0</b>	2	45	Blue	12		8886442043	
	<b>3-0</b>	2	75	Blue	36		88864554-41	
<b>C-14</b> 3/8 Circle 24 mm 	<b>4-0</b>	1.5	45	Blue	12		8886442233	
	<b>3-0</b>	2	45	Blue	12		8886442243	
	<b>3-0</b>	2	90	Blue	36		8886400241	
	<b>2-0</b>	3	45	Blue	12		8886442253	
	<b>2-0</b>	3	75	Blue	12		8886442153	






## REVERSE CUTTING

Needles Specifications	Size		Length cm	Color	Units Per box	Comments	Reorder code
	USP	EP					
<b>C-16</b> 3/8 Circle    30 mm 	<b>4-0</b>	1.5	75	Blue	36		8886442431
	<b>3-0</b>	2	75	Blue	36		8886442441
	<b>2-0</b>	3	75	Blue	36		8886442451
	<b>0</b>	3.5	75	Blue	36		8886442461
<b>C-17</b> 3/8 Circle    39 mm 	<b>3-0</b>	2	75	Blue	36		8886442641
	<b>2-0</b>	3	75	Blue	36		8886442651
	<b>0</b>	3.5	75	Blue	36		8886442661
<b>GS-11</b> 1/2 Circle    37 mm 	<b>0</b>	3.5	100	Blue	12		8886400963
	<b>1</b>	4	75	Blue	36		8886443371
<b>GS-13</b> 1/2 Circle    48 mm (Not to scale) 	<b>1</b>	4	75	Blue	36		8886443671
	<b>1</b>	4	100	Blue	12		CPB-738G

## PREMIUM REVERSE CUTTING

Needles Specifications	Size USP	EP	Length cm	Color	Units Per box	Comments	Reorder code
<b>P-10</b> 3/8 Circle 11 mm 	<b>7-0</b>	0.5	45	Blue	12	Surgalloy™	8886441003
	<b>6-0</b>	0.7	45	Blue	12	Surgalloy™	8886441013
	<b>5-0</b>	1	45	Blue	12	Surgalloy™	8886441023
<b>P-13</b> 3/8 Circle 13 mm 	<b>6-0</b>	0.7	45	Blue	12	Surgalloy™	SPB-1213G
	<b>5-0</b>	1	45	Blue	12	Surgalloy™	SPB-5223G
	<b>5-0</b>	1	75	Blue	36	Surgalloy™	8886402721
	<b>4-0</b>	1.5	45	Blue	12	Surgalloy™	SPB-1233G
<b>SBE-2</b> 3/8 Circle 13 mm 	<b>6-0</b>	0.7	45	Blue	12		8886440013
	<b>5-0</b>	1	45	Blue	12		8886440023
<b>P-11</b> 3/8 Circle 16 mm 	<b>6-0</b>	0.7	45	Blue	12	Surgalloy™	SPB-5413G
	<b>5-0</b>	1	45	Blue	12	Surgalloy™	SPB-5142G
	<b>5-0</b>	1	75	Blue	36	Surgalloy™	8886402821
	<b>4-0</b>	1.5	45	Blue	12	Surgalloy™	SPB-5143G
	<b>4-0</b>	1.5	45	Blue	36	Surgalloy™	8886441431
	<b>4-0</b>	1.5	75	Blue	36	Surgalloy™	8886402831
<b>SBE-3</b> 3/8 Circle 16 mm 	<b>6-0</b>	0.7	45	Blue	12		8886440113
	<b>5-0</b>	1	45	Blue	12		8886440123
	<b>4-0</b>	1.5	45	Blue	12		8886440133

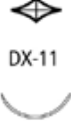
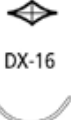

## PREMIUM REVERSE CUTTING

Needles Specifications	Size USP	EP	Length cm	Color	Units Per box	Comments	Reorder code
<b>P-12</b> 3/8 Circle 19 mm 	5-0	1	45	Blue	12	Surgalloy™	SPB-1623G
	4-0	1.5	45	Blue	12	Surgalloy™	SPB-5633G
	4-0	1.5	75	Blue	36	Surgalloy™	8886442531
	3-0	2	45	Blue	12	Surgalloy™	SPB-1643G
	3-0	2	45	Blue	36	Surgalloy™	8886441543
<b>SBE-4</b> 3/8 Circle 19 mm 	5-0	1	45	Blue	12		8886440223
	4-0	1.5	45	Blue	12		8886440233
	3-0	2	45	Blue	12		8886440243
<b>P-14</b> 3/8 Circle 24 mm 	4-0	1.5	45	Blue	12	Surgalloy™	SPB-5833G
	3-0	2	45	Blue	12	Surgalloy™	SPB-5843G
	3-0	2	75	Blue	36	Surgalloy™	88864582-43
	2-0	3	45	Blue	12	Surgalloy™	88864418-53
	2-0	3	75	Blue	36	Surgalloy™	88864582-53
<b>SBE-6</b> 3/8 Circle 24 mm 	4-0	1.5	45	Blue	12		8886440333
	3-0	2	45	Blue	12		8886440343
<b>P-24</b> 1/2 Circle 16 mm 	4-0	1.5	45	Clear	12	Surgalloy™	SPB-5433G







## CONVENTIONAL CUTTING

Needles Specifications	Size		Length cm	Color	Units Per box	Comments	Reorder code
	USP	EP					
<b>PC-12</b> 3/8 Circle    19 mm 	<b>5-0</b>	1	45	Blue	12		SPB-1945G
<b>SC-2</b> Straight    60 mm 	<b>3-0</b>	2	75	Blue	36		88864490-41
	<b>2-0</b>	3	75	Blue	36		88864490-51

## X CUTTING

Needles Specifications	Size		Length cm	Color	Units Per box	Comments	Reorder code
	USP	EP					
<b>DX-11</b> 3/8 Circle    11 mm 	<b>6-0</b>	0.7	45	Blue	12	DermaX	XNF-697
<b>DX-16</b> 3/8 Circle    16 mm 	<b>6-0</b>	0.7	45	Blue	12	DermaX	XNF-1966
	<b>5-0</b>	1	45	Blue	12	DermaX	XNF-1965
	<b>4-0</b>	1.5	45	Blue	12	DermaX	XNF-1964
<b>DX-19</b> 3/8 Circle    19 mm 	<b>5-0</b>	1	45	Blue	12	DermaX	XNF-1995
	<b>4-0</b>	1.5	45	Blue	12	DermaX	XNF-1994
	<b>3-0</b>	2	45	Blue	12	DermaX	XNF-1993

## TAPER POINT


Needles Specifications	Size USP	EP	Length cm	Color	Units Per box	Comments	Reorder code
<b>V-20</b> 1/2 Circle 26 mm 	<b>3-0</b>	2	75	Blue	36	Surgalloy™	8886445241
	<b>2-0</b>	3	75	Blue	36	Surgalloy™	8886445251
	<b>0</b>	3.5	75	Blue	36	Surgalloy™	8886445261
<b>GS-22</b> 1/2 Circle 27 mm 	<b>2-0</b>	3	75	Blue	36		8886445951
	<b>2-0</b>	3	5×45	Blue	12	D-Tach™	8886445853
	<b>0</b>	3.5	5×45	Blue	12	D-Tach™	8886445863
	<b>0</b>	3.5	75	Blue	36		8886445961
	<b>1</b>	4	75	Blue	36		8886445971
<b>GS-22 LOOP</b> 1/2 Circle 27 mm 	<b>2-0</b>	3	150	Blue	36	Loop	88864619-51
<b>V-26</b> 1/2 Circle 37 mm 	<b>3-0</b>	2	75	Blue	36		8886445041
<b>GS-21</b> 1/2 Circle 37 mm 	<b>0</b>	3.5	5×45	Blue	12	D-Tach™	8886445463
	<b>0</b>	3.5	75	Blue	36		8886445561
	<b>1</b>	4	5×45	Blue	12	D-Tach™	8886445473
	<b>1</b>	4	75	Blue	36		8886445571
<b>GS-24</b> 1/2 Circle 40 mm 	<b>1</b>	4	75	Blue	36		88864463-71

## TAPER POINT

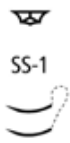

Needles Specifications	Size		Length		Color	Units Per box	Comments	Reorder code
	USP	EP	cm	Inches				
<b>GS-25</b> 1/2 Circle 48 mm (Not to scale) 	<b>2</b>	5	75	30	Blue	36		88864465-81
	<b>1</b>	4	5×45	18	Blue	12	D-Tach™	88864464-73
	<b>1</b>	4	75	30	Blue	36		88864465-71
<b>GS-25 LOOP</b> 1/2 Circle 48 mm 	<b>0</b>	3.5	200	80	Blue	12	Loop	88864538-63
<b>GS-26</b> 1/2 Circle 65 mm 	<b>2</b>	5	150	60	Blue	36		88864470-81
<b>GS-26 LOOP</b> 1/2 Circle 65 mm 	<b>0</b>	3.5	150	60	Blue	36	Loop	8886440861



## HEAVY TAPER POINT

Needles Specifications	Size		Length		Color	Units Per box	Comments	Reorder code
	USP	EP	cm	Inches				
<b>HGS-22</b> 1/2 Circle 27 mm 	<b>0</b>	3.5	5×45	18	Blue	12	D-Tach™	8886446063

## SPATULA

Needles Specifications	Size		Length		Color	Units Per box	Comments	Reorder code
	USP	EP	cm	Inches				
<b>SS-1</b> 1/6 Circle 5.6 mm Double-arm Curvature 60° 	<b>5-0</b>	1	30	12	Blue	12	Surgalloy™	PB-7740K
<b>SS-24</b> 1/4 Circle 8.6 mm Double-arm Curvature 95° 	<b>6-0</b>	0.7	30	12	Blue	12	Surgalloy™	PB-6713K
	<b>5-0</b>	1	30	12	Blue	12	Surgalloy™	PB-6723K

# SURGILON™

Braided Nylon

Coated with Silicone to Enhance Handling.



Please refer to the individual package IFUs for complete instructions, indications, contraindications, warnings and precautions.

<b>Description</b>	SURGILON™ braided sutures are inert, nonabsorbable, sterile sutures composed of the long-chain aliphatic polymers Nylon 6 and Nylon 6.6. The braided sutures are coated uniformly with silicone to enhance handling characteristics, ease of passage through tissue and reduction of capillarity.
<b>Indications</b>	SURGILON™ braided nylon sutures are indicated for use in general soft tissue approximation and/or ligation, including use in cardiovascular, ophthalmic, and neurological surgery.
<b>Actions</b>	SURGILON™ braided nylon sutures all elicit a minimal acute inflammatory reaction in tissue, which is followed by gradual encapsulation of the suture by fibrous connective tissue.
<b>Tensile Strength</b>	As with any nylon suture, though not absorbed, progressive hydrolysis of the suture, in vivo, may result in gradual loss of its tensile strength over time.
<b>Sterilization</b>	Surgilon™ braided nylon sutures are sterilized by ethylene oxide in sizes 6-0 (0.7 metric) and smaller and by gamma irradiation in sizes 2 (5 metric) through 5-0 (1 metric).
<b>Packaging</b>	SURGILON™ braided nylon sutures are available dyed black in sizes 1 (4 Metric) through 4-0 (1.5 Metric) and size 6-0 (0.7 metric). The sutures are supplied sterile, in pre-cut lengths or affixed to needles using both permanent and removable needle attachment techniques.






**PRE-CUT**

Needles Specifications	Size		Length cm	Color	Units Per box	Comments	Reorder code
	USP	EP					
	<b>4-0</b>	1.5	12×45	Black	36	Pre-cut, RACETRAC	88861918-31
	<b>4-0</b>	1.5	7×75	Black	36	Pre-cut	88861919-31
	<b>3-0</b>	2	12×45	White	36	Pre-cut, RACETRAC	88861899-41
	<b>3-0</b>	2	12×45	Black	36	Pre-cut, RACETRAC	88861918-41
	<b>3-0</b>	2	7×75	Black	36	Pre-cut	88861919-41
	<b>2-0</b>	3	12×45	White	36	Pre-cut, RACETRAC	88861899-51
	<b>2-0</b>	3	12×45	Black	36	Pre-cut, RACETRAC	88861918-51
	<b>2-0</b>	3	7×75	Black	36	Pre-cut	88861919-51
	<b>0</b>	3.5	6×45	Black	36	Pre-cut, RACETRAC	88861917-61
	<b>0</b>	3.5	6×45	White	36	Pre-cut, RACETRAC	88861902-61
	<b>0</b>	3.5	7×75	Black	36	Pre-cut	88861919-61
	<b>1</b>	4	6×45	Black	36	Pre-cut, RACETRAC	88861917-71
	<b>1</b>	4	6×45	White	36	Pre-cut, RACETRAC	88861902-71
	<b>1</b>	4	7×75	Black	36	Pre-cut	88861919-71
	<b>2</b>	5	6×45	Black	24	Pre-cut, RACETRAC	88861917-82



**REEL**

Needles Specifications	Size		Length cm	Color	Units Per box	Comments	Reorder code
	USP	EP					
	<b>3-0</b>	2	366	Black	24	Reel	88861906-42
	<b>2-0</b>	3	366	Black	24	Reel	88861906-52
	<b>0</b>	3.5	250	Black	24	Reel	88861920-62
	<b>1</b>	4	250	Black	24	Reel	88861920-72
	<b>2</b>	5	300	Black	12	Reel	88861914-83
	<b>3</b>	6	238	Black	12	Reel	88861920-93




## REVERSE CUTTING

Needles Specifications			Size		Length cm	Color	Units Per box	Comments	Reorder code
			USP	EP					
<b>C-13</b> 3/8 Circle 19 mm		<b>3-0</b>	2	45	Black	12		88861931-43	
		<b>2-0</b>	3	75	Black	36		88861982-51	
<b>C-14</b> 3/8 Circle 24 mm		<b>3-0</b>	2	75	Black	36		88861935-41	
		<b>2-0</b>	3	45	Black	12		88861934-53	
		<b>2-0</b>	3	75	Black	36		88861935-51	
		<b>0</b>	3.5	75	Black	36		88861935-61	
<b>C-16</b> 3/8 Circle 30 mm		<b>3-0</b>	2	5 × 45	Black	24	D-Tach™	88861991-42	
		<b>2-0</b>	3	75	Black	36		88861883-51	
		<b>0</b>	3.5	75	Black	36		88861883-61	
		<b>1</b>	4	75	Black	36		88861883-71	
<b>C-18</b> 3/8 Circle 35 mm		<b>1</b>	4	75	Black	36		88861922-71	
<b>GS-10</b> 1/2 Circle 26 mm		<b>2-0</b>	3	75	Black	36		88861932-51	



## REVERSE CUTTING

Needles Specifications	Size		Length cm	Color	Units Per box	Comments	Reorder code
	USP	EP					
<b>HOS-10</b> 1/2 Circle 26 mm 	1	4	75	Black	36		88861896-71
<b>GS-11</b> 1/2 Circle 37 mm 	0	3.5	75	Black	36		88861915-61
	1	4	75	Black	36		88861915-71
<b>HOS-12</b> 1/2 Circle 40 mm 	1	4	5 × 45	Black	24	D-Tach™	88861913-72
	1	4	75	Black	36		88861895-71




## PREMIUM REVERSE CUTTING

Needles Specifications	Size		Length cm	Color	Units Per box	Comments	Reorder code
	USP	EP					
<b>P-13</b> 3/8 Circle 13 mm 	5-0	1	45	Black	12	Surgalloy™	SBS-1928G
	4-0	1.5	45	White	12	Surgalloy™	SBS-1880G
<b>P-11</b> 3/8 Circle 16 mm 	4-0	1.5	45	Black	12		88861981-33
<b>P-12</b> 3/8 Circle 19 mm 	4-0	1.5	45	White	12	Surgalloy™	SBS-1884G
	3-0	2	45	White	12	Surgalloy™	SBS-1883G

## CONVENTIONAL CUTTING

Needles Specifications	Size		Length cm	Color	Units Per box	Comments	Reorder code
	USP	EP					
<b>SC-1</b> Straight      51 mm      ▲ 	<b>2-0</b>	3	75	Black	36		88861985-51
	<b>0</b>	3.5	75	Black	36		88861985-61
	<b>1</b>	4	75	Black	36		88861985-71
<b>SC-2</b> Straight      60 mm      ▲ 	<b>0</b>	3.5	75	Black	36		88861948-61
	<b>2</b>	5	75	Black	36		88861948-81

## TAPER POINT




Needles Specifications	Size		Length cm	Color	Units Per box	Comments	Reorder code
	USP	EP					
<b>CV-20</b> 1/2 Circle      10 mm      ⊙ 	<b>5-0</b>	1	5 × 45	Black	24	D-Tach™	88861910-22
	<b>4-0</b>	1.5	5 × 45	Black	24	D-Tach™	88861910-32
<b>CV-22</b> 1/2 Circle      13 mm      ⊙ 	<b>4-0</b>	1.5	5 × 45	Black	24	D-Tach™	88861957-32
<b>CV-23</b> 1/2 Circle      17 mm      ⊙ 	<b>4-0</b>	1.5	5 × 45	Black	24	D-Tach™	88861927-32
	<b>4-0</b>	1.5	75	Black	36		88861995-31
	<b>3-0</b>	2	5 × 45	Black	24	D-Tach™	88861927-42
	<b>3-0</b>	2	75	Black	36		88861995-41

## TAPER POINT


Needles Specifications	Size USP	EP	Length cm	Color	Units Per box	Comments	Reorder code
<b>CV-24</b> 1/2 Circle 20 mm 	4-0	1.5	5 × 45	Black	24	D-Tach™	88861989-32
	3-0	2	5 × 45	Black	24	D-Tach™	88861989-42
<b>GS-23</b> 1/2 Circle 22 mm 	0	3.5	5 × 45	Black	24	D-Tach™	88861993-62
<b>V-20</b> 1/2 Circle 26 mm 	3-0	2	5 × 45	Black	24	D-Tach™	88861958-42
	3-0	2	75	Black	36		88861951-41
	3-0	2	75	Black	24	D-Tach™	88861926-42
	2-0	3	45	Black	24	D-Tach™	88861958-52
	2-0	3	5 × 75	Black	24	D-Tach™	88861967-52
	2-0	3	75	Black	36		88861951-51
	0	3.5	75	Black	36		88861951-61
<b>GS-22</b> 1/2 Circle 27 mm 	2-0	3	5 × 45	Black	24	D-Tach™	88861969-52
	0	3.5	5 × 45	Black	24	D-Tach™	88861969-62
	0	3.5	75	Black	36		88861968-61
	1	4	5 × 45	Black	24	D-Tach™	88861969-72
<b>V-30</b> 1/2 Circle 30 mm 	2-0	3	75	Black	36		88861964-51
	0	3.5	75	Black	36		88861964-61




## TAPER POINT

Needles Specifications	Size		Length cm	Color	Units Per box	Comments	Reorder code
	USP	EP					
<b>GS-21</b> 1/2 Circle 37 mm 	<b>2-0</b>	3	5 × 45	Black	24	D-Tach™	88861962-52
	<b>0</b>	3.5	75	Black	36		88861971-61
	<b>1</b>	4	5 × 45	Black	24	D-Tach™	88861962-72
	<b>1</b>	4	75	Black	36		88861971-71
<b>GS-25</b> 1/2 Circle 48 mm 	<b>0</b>	3.5	5 × 45	Black	24	D-Tach™	88861998-62
<b>SK</b> SKI 23 mm 	<b>2-0</b>	3	250	Black	36		88862011-51

## HEAVY TAPER POINT

Needles Specifications	Size		Length cm	Color	Units Per box	Comments	Reorder code
	USP	EP					
<b>HGS-22</b> 1/2 Circle 27 mm 	<b>2-0</b>	3	5 × 45	Black	24	D-Tach™	88861928-52
	<b>0</b>	3.5	5 × 45	Black	24	D-Tach™	88861928-62
	<b>1</b>	4	5 × 45	Black	24	D-Tach™	88861928-72

## DIAMOND POINT

Needles Specifications	Size		Length cm	Color	Units Per box	Comments	Reorder code
	USP	EP					
<b>KV-25</b> 1/2 Circle 22 mm 	<b>0</b>	3.5	5 × 45	Black	24	D-Tach™	88861909-62

# SOFSILK™

Wax Coated Braided Silk

Maximum Surgical Handling.



Please refer to the individual package IFUs for complete instructions, indications, contraindications, warnings and precautions.

**Description**

Sofsilk™ braided silk sutures are nonabsorbable, sterile, non-mutagenic surgical sutures composed of natural proteinaceous silk fibers called fibroin. This protein is derived from the domesticated silkworm species *Bombyx mori* of the family *bombycidae*. The silk fibers are treated to remove the naturally- occurring sericin gum and braided to produce Sofsilk™ surgical silk sutures. The braided sutures are available coated uniformly with either silicone or a special wax mixture to reduce capillarity and to increase surface lubricity which enhances handling characteristics, ease of passage through tissue, and knot run-down properties.

**Indications**

Sofsilk™ sutures are indicated for use in general soft tissue approximation and/or ligation, including use in cardiovascular, ophthalmic, microsurgery, and neurological surgery.

**Actions**

Sofsilk™ sutures elicit a minimal acute inflammatory reaction in tissue, which is followed by gradual encapsulation of the suture by fibrous connective tissue.

**Tensile Strength**

While Sofsilk™ sutures are not absorbed, progressive degradation of the proteinaceous silk fiber in vivo may result in gradual loss of the suture's tensile strength over time.

**Sterilization**

Gamma irradiation.

**Packaging**

Sofsilk™ sutures are available in black with Logwood extract and uniformly coated with silicone or a special wax mixture, in sizes 2 (5 Metric) through 8-0 (0.4 Metric). The sutures are supplied sterile, in pre-cut lengths and ligating reels, nonneedled or affixed to needles using both permanent and removable needle attachment techniques. Please refer to the individual package IFUs for complete instructions, indications, contraindications, warnings and precautions.

PRE-CUT

Needles Specifications	Size		Length cm	Color	Units Per box	Comments	Reorder code
	USP	EP					
	5-0	1	12×45	Black	36	Pre-cut	S-182
	4-0	1.5	12×45	Black	24	Pre-cut, racetrac	S-193
	4-0	1.5	12×45	Black	36	Pre-cut	S-183
	4-0	1.5	12×60	Black	36	Pre-cut	S-243
	4-0	1.5	12×75	Black	36	Pre-cut	S-303
	3-0	2	12×45	Black	24	Pre-cut, racetrac	S-194
	3-0	2	12×45	Black	36	Pre-cut	S-184
	3-0	2	12×60	Black	24	Pre-cut, racetrac	S-254
	3-0	2	12×60	Black	36	Pre-cut	S-244
	3-0	2	12×75	Black	24	Pre-cut, racetrac	S-404
	3-0	2	12×75	Black	36	Pre-cut	S-304
	2-0	3	12×45	Black	24	Pre-cut, racetrac	S-195
	2-0	3	12×45	Black	36	Pre-cut	S-185
	2-0	3	12×60	Black	36	Pre-cut	S-245
	2-0	3	12×75	Black	36	Pre-cut	S-305
	2-0	3	2×150	Black	36	Pre-cut	S-605
	0	3.5	6×45	Black	24	Pre-cut, racetrac	S-196
	0	3.5	12×45	Black	24	Pre-cut	S-176
	0	3.5	12×60	Black	24	Pre-cut	S-346
	0	3.5	6×60	Black	36	Pre-cut	S-246
	0	3.5	12×75	Black	24	Pre-cut	S-316
	0	3.5	2×150	Black	36	Pre-cut	S-606
	1	4	6×45	Black	24	Pre-cut, racetrac	S-197
	1	4	6×45	Black	36	Pre-cut	S-187
	1	4	12×75	Black	24	Pre-cut	S-317
	1	4	2×150	Black	36	Pre-cut	S-607
	2	5	12×75	Black	24	Pre-cut	S-318
	2	5	2×150	Black	36	Pre-cut	S-608

## PRE-CUT

Needles Specifications	Size		Length cm	Color	Units Per box	Comments	Reorder code
	USP	EP					
	5	7	3 × 60	Black	36	Pre-cut	S-282
	5	7	2 × 150	Black	12	Pre-cut	S-610

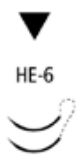

## NON-NEEDED-STANDARD LENGTH

Needles Specifications	Size		Length cm	Color	Units Per box	Comments	Reorder code
	USP	EP					
	2	5	100	Black	36	Standard length	S-206
	4	6	180	Black	12	Standard length	S-199




## REEL

Needles Specifications	Size		Length cm	Color	Units Per box	Comments	Reorder code
	USP	EP					
	4-0	1.5	250	Black	24	Reel	LS-640
	3-0	2	250	Black	24	Reel	LS-639
	2-0	3	250	Black	24	Reel	LS-638
	0	3.5	250	Black	24	Reel	LS-637
	1	4	250	Black	24	Reel	LS-636



## REVERSE CUTTING

Needles Specifications	Size		Length cm	Color	Units Per box	Comments	Reorder code	
	USP	EP						
<b>HE-6</b> 3/8 Circle 7.4 mm Double-arm Curvature 110°		7-0	0.5	45	Black	12	S-1768K	
		6-0	0.7	45	Black	12	Surgalloy™	S-1769K
<b>HE-5</b> 3/8 Circle 8.6 mm Double-arm Curvature 135°		6-0	0.7	30	Black	12	Surgalloy™	S-1746K




## REVERSE CUTTING

Needles Specifications	Size USP	EP	Length cm	Color	Units Per box	Comments	Reorder code	
<b>HE-1</b> 3/8 Circle 11 mm Double-arm Curvature 135°	6-0	0.7	45	Black	12	Surgalloy™	S-1750K	
<b>C-1</b> 3/8 Circle 12 mm		6-0	0.7	45	Black	36	S-1172	
		5-0	1	45	Black	36	S-1173	
		4-0	1.5	45	Black	36	S-1174	
		4-0	1.5	75	Black	36	S-1373	
<b>C-1</b> 3/8 Circle 12 mm Double-arm		6-0	0.7	45	Black	12	S-1272	
		4-0	1.5	30	Black	12	S-1735	
		4-0	1.5	45	Black	12	S-1274	
<b>HE-3</b> 3/8 Circle 13 mm Double-arm Curvature 155°		4-0	1.5	45	Black	12	Surgalloy™	S-1789K
		4-0	1.5	45	Black	12	Surgalloy™	S-1783K

## REVERSE CUTTING


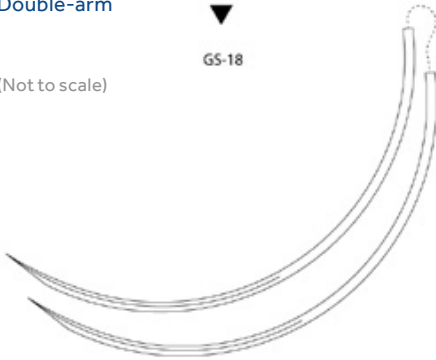


Needles Specifications	Size USP	EP	Length cm	Color	Units Per box	Comments	Reorder code
<b>C-12</b> 3/8 Circle 16 mm 	6-0	0.7	45	Black	36		SS-681
	5-0	1	75	Black	36		SS-673
	4-0	1.5	45	Black	12		SS-686
	4-0	1.5	75	Black	36		SS-645
	3-0	2	75	Black	36		SS-745
	2-0	3	75	Black	36		SS-746
<b>C-13</b> 3/8 Circle 19 mm 	5-0	1	45	Black	36		SS-682
	5-0	1	75	Black	36		SS-652
	4-0	1.5	45	Black	12		SS-683G
	4-0	1.5	45	Black	36		SS-683
	4-0	1.5	75	Black	36		SS-653
	3-0	2	45	Black	36		SS-694
	3-0	2	75	Black	36		SS-654
	2-0	3	75	Black	36		SS-655
	0	3.5	75	Black	36		SS-656

## REVERSE CUTTING



Needles Specifications	Size		Length cm	Color	Units Per box	Comments	Reorder code
	USP	EP					
<b>C-14</b> 3/8 Circle    24 mm 	4-0	1.5	45	Black	36		SS-629
	3-0	2	45	Black	36		SS-684
	3-0	2	45	Black	12		SS-684G
	3-0	2	75	Black	36		SS-784
	2-0	3	45	Black	36		SS-695
	2-0	3	75	Black	36		SS-785
	0	3.5	75	Black	36		SS-786
	1	4	75	Black	36		SS-787
<b>C-15</b> 3/8 Circle    26 mm 	2-0	3	45	Black	12		SS-685G
	2-0	3	45	Black	36		SS-685
<b>C-16</b> 3/8 Circle    30 mm 	3-0	2	75	Black	36		SS-675
	2-0	3	45	Black	36		SS-677
	2-0	3	75	Black	36		SS-679
	0	3.5	45	Black	36		SS-678
	0	3.5	5 × 75	Black	24	D-Tach™	SS-70-M
	0	3.5	75	Black	36		SS-680
	1	4	75	Black	36		SS-689









## REVERSE CUTTING

Needles Specifications	Size		Length cm	Color	Units Per box	Comments	Reorder code
	USP	EP					
<b>C-17</b> 3/8 Circle    39 mm 	<b>2-0</b>	3	75	Black	36		SS-648
	<b>0</b>	3.5	75	Black	36		SS-647
	<b>1</b>	4	75	Black	36		SS-646
	<b>2</b>	5	75	Black	36		SS-651
<b>GS-18</b> 3/8 Circle    77 mm Double-arm (Not to scale) 	<b>2</b>	5	75	Black	24		CS-482
	<b>5</b>	7	75	Black	24		CS-485
<b>HE-7</b> 1/2 Circle    7.4 mm Double-arm Curvature 145° 	<b>7-0</b>	0.5	45	Black	12		S-1766K
	<b>6-0</b>	0.7	45	Black	12	Surgalloy™	S-1765K
<b>C-21</b> 1/2 Circle    13 mm 	<b>4-0</b>	1.5	5 × 45	Black	24	D-Tach™	SS-82-M

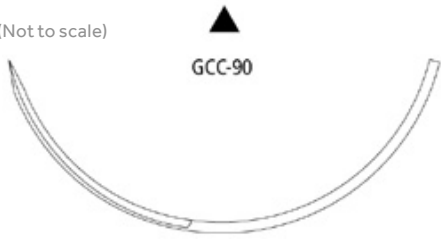

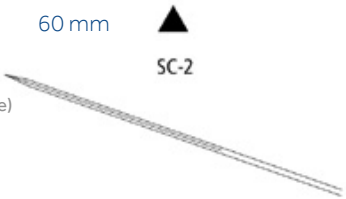
## REVERSE CUTTING

Needles Specifications	Image	Size		Length cm	Color	Units Per box	Comments	Reorder code
		USP	EP					
<b>C-22</b> 1/2 Circle    20 mm		4-0	1.5	45	Black	36		SS-733
<b>C-23</b> 1/2 Circle    24 mm		3-0	2	45	Black	36		SS-732
		3-0	2	75	Black	36		SS-632
		2-0	3	45	Black	36		SS-734
		2-0	3	75	Black	36		SS-633
<b>GS-10</b> 1/2 Circle    26 mm		2-0	3	45	Black	36		CS-390
<b>GS-11</b> 1/2 Circle    37 mm		2-0	3	75	Black	36		CS-490

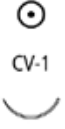
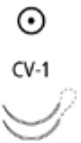
## PREMIUM REVERSE CUTTING

Needles Specifications	Size USP	EP	Length cm	Color	Units Per box	Comments	Reorder code
<b>P-10</b> 3/8 Circle 11 mm 	6-0	0.7	45	Black	36	Surgalloy™	SS-5639
	6-0	0.7	45	Black	12	Surgalloy™	SS-1639G
<b>P-13</b> 3/8 Circle 13 mm 	6-0	0.7	45	Black	12	Surgalloy™	SS-5649G
	5-0	1	45	Black	36	Surgalloy™	SS-5640
	4-0	1.5	45	Black	12	Surgalloy™	SS-5641G
	4-0	1.5	45	Black	36	Surgalloy™	SS-5641
<b>P-11</b> 3/8 Circle 16 mm 	5-0	1	45	Black	36	Surgalloy™	SS-5678
<b>P-12</b> 3/8 Circle 19 mm 	5-0	1	45	Black	36	Surgalloy™	SS-5676
	4-0	1.5	45	Black	36	Surgalloy™	SS-5677
	3-0	2	45	Black	36	Surgalloy™	SS-5679
<b>P-14</b> 3/8 Circle 24 mm 	3-0	2	45	Black	36	Surgalloy™	SS-5684
	2-0	3	45	Black	12		SS-5685G
<b>P-24</b> 1/2 Circle 16 mm 	4-0	1.5	45	Black	12	Surgalloy™	SS-1694G

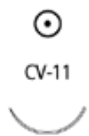


## CONVENTIONAL CUTTING

Needles Specifications	Size		Length cm	Color	Units Per box	Comments	Reorder code
	USP	EP					
<b>GCC-90</b> 3/8 Circle 90 mm (Not to scale) 	<b>2-0</b>	3	100	Black	36		CS-791
	<b>0</b>	3.5	100	Black	36		CS-792
	<b>1</b>	4	100	Black	36		CS-793
	<b>2</b>	5	100	Black	36		CS-794
<b>SC-1</b> Straight 51 mm Double-arm (Not to scale) 	<b>2-0</b>	3	75	Black	36		SS-523
	<b>1</b>	4	75	Black	36		SS-525
<b>SC-2</b> Straight 60 mm (Not to scale) 	<b>4-0</b>	1.5	75	Black	36		SS-621
	<b>3-0</b>	2	75	Black	36		SS-622
	<b>3-0</b>	2	100	Black	36		SS-722
	<b>2-0</b>	3	75	Black	36		SS-623
	<b>2-0</b>	3	100	Black	36		SS-723
	<b>0</b>	3.5	75	Black	36		SS-624





## TAPER POINT

Needles Specifications	Size		Length cm	Color	Units Per box	Comments	Reorder code
	USP	EP					
<b>CV-1</b> 3/8 Circle 9 mm 	<b>6-0</b>	0.7	45	Black	36	Surgalloy™	VS-802
	<b>7-0</b>	0.5	45	Black	36	Surgalloy™	VS-809
<b>CV-1</b> 3/8 Circle 9 mm Double-arm 							

## TAPER POINT

Needles Specifications	Size USP	EP	Length cm	Color	Units Per box	Comments	Reorder code
<b>CV-11</b> 3/8 Circle 13 mm 	6-0	0.7	75	Black	36	Surgalloy™	VS-889
	5-0	1	75	Black	36	Surgalloy™	VS-890
	4-0	1.5	75	Black	36	Surgalloy™	VS-891
<b>CV-11</b> 3/8 Circle 13 mm Double-arm 	6-0	0.7	45	Black	36	Surgalloy™	VS-806
	6-0	0.7	5 × 60	Black	12	Surgalloy™	VS-706-M
	4-0	1.5	45	Black	36	Surgalloy™	VS-709
<b>CV-15</b> 3/8 Circle 17 mm 	5-0	1	75	Black	36	Surgalloy™	VS-880
	4-0	1.5	75	Black	36	Surgalloy™	VS-881
	3-0	2	75	Black	36	Surgalloy™	VS-882
	2-0	3	75	Black	36	Surgalloy™	VS-823
<b>CV-15</b> 3/8 Circle 17 mm Double-arm 	4-0	1.5	60	Black	36	Surgalloy™	VS-581




## TAPER POINT

Needles Specifications	Size USP	EP	Length cm	Color	Units Per box	Comments	Reorder code
<b>CVF-21</b> 1/2 Circle 12 mm  CVF-21 	5-0	1	45	Black	36	Surgalloy™	CS-210
	4-0	1.5	45	Black	36	Surgalloy™	CS-211
<b>CV-22</b> 1/2 Circle 13 mm  CV-22 	4-0	1.5	5 × 45	Black	24	D-Tach™	GS-44-M

## TAPER POINT

Needles Specifications	Size USP	EP	Length cm	Color	Units Per box	Comments	Reorder code
<b>CV-23</b> 1/2 Circle 17 mm 	5-0	1	75	Black	36	Surgalloy™	VS-870
	4-0	1.5	75	Black	36	Surgalloy™	VS-871
	3-0	2	5 × 45	Black	24	D-Tach™	GS-47-M
	3-0	2	75	Black	36	Surgalloy™	VS-872
	2-0	3	75	Black	36	Surgalloy™	VS-873
<b>CV-24</b> 1/2 Circle 20 mm 	4-0	1.5	75	Black	36		GS-451
	3-0	2	75	Black	36		GS-452
	2-0	3	75	Black	36		GS-453
<b>CV-25</b> 1/2 Circle 22 mm 	4-0	1.5	5 × 45	Black	24	D-Tach™	GS-34-M
	3-0	2	5 × 45	Black	24	D-Tach™	GS-33-M
	3-0	2	5 × 75	Black	24	D-Tach™	GS-30-M
<b>CV-25</b> 1/2 Circle 22 mm Double-arm 	3-0	2	75	Black	36	Surgalloy™	VS-552
<b>GS-23</b> 1/2 Circle 22 mm 	2-0	3	75	Black	36		CS-562

## TAPER POINT




Needles Specifications	Size		Length cm	Color	Units Per box	Comments	Reorder code
	USP	EP					
<b>V-20</b> 1/2 Circle 26 mm 	4-0	1.5	5 × 45	Black	24	D-Tach™	GS-64-M
	4-0	1.5	75	Black	36		GS-831
	3-0	2	5 × 45	Black	24	D-Tach™	GS-63-M
	3-0	2	5 × 75	Black	24	D-Tach™	GS-67-M
	3-0	2	75	Black	36		GS-832
	2-0	3	5 × 45	Black	24	D-Tach™	GS-62-M
	2-0	3	75	Black	36		GS-833
	2-0	3	5 × 75	Black	24	D-Tach™	GS-66-M
	0	3.5	75	Black	36		GS-834
	1	4	75	Black	36		GS-835
<b>V-20</b> 1/2 Circle 26 mm Double-arm 	2-0	3	90	Black	36	Surgalloy™	VS-533
<b>GS-22</b> 1/2 Circle 27 mm 	2-0	3	5 × 45	Black	24	D-Tach™	CS-16-M
	0	3.5	5 × 45	Black	24	D-Tach™	CS-17-M



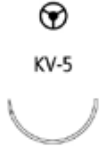


## TAPER POINT

Needles Specifications	Size		Length cm	Color	Units Per box	Comments	Reorder code
	USP	EP					
<b>V-30</b> 1/2 Circle 30 mm 	<b>3-0</b>	2	75	Black	36		GS-822
	<b>2-0</b>	3	75	Black	36		GS-823
	<b>0</b>	3.5	75	Black	36		GS-824
<b>GS-21</b> 1/2 Circle 37 mm 	<b>2-0</b>	3	75	Black	36		CS-423
	<b>0</b>	3.5	5 × 45	Black	24	D-Tach™	CS-10-M
	<b>0</b>	3.5	75	Black	36		CS-424
	<b>1</b>	4	75	Black	36		CS-425
<b>V-26</b> 1/2 Circle 37 mm 	<b>3-0</b>	2	75	Black	36		VS-842
	<b>2-0</b>	3	75	Black	36		VS-843
	<b>0</b>	3.5	75	Black	36		VS-844
	<b>1</b>	4	75	Black	36		VS-845
	<b>1</b>	4	2 × 75	White	24	D-Tach™	CS-745
	<b>1</b>	4	75	White	36		CS-744
	<b>2</b>	5	75	Black	36		VS-846
	<b>2</b>	5	2 × 75	White	24	D-Tach™	CS-746





## HEAVY TAPER POINT

Needles Specifications	Size USP	EP	Length cm	Color	Units Per box	Comments	Reorder code
<b>HGS-23</b> 1/2 Circle 22 mm 	2-0	3	5 × 45	Black	24	D-Tach™	CS-93-M
	0	3.5	5 × 45	Black	24	D-Tach™	CS-92-M
<b>HGS-22</b> 1/2 Circle 27 mm 	0	3.5	5 × 45	Black	24	D-Tach™	CS-85-M
<b>HGS-20</b> 1/2 Circle 34mm 	1	4	5 × 75	Black	12		X-1095

## TAPERCUTTING

Needles Specifications	Size USP	EP	Length cm	Color	Units Per box	Comments	Reorder code
<b>KV-5</b> 1/2 Circle 17 mm 	3-0	2	75	Black	36	Surgalloy™	VS-863
<b>KV-7</b> 1/2 Circle 26 mm 	2-0	3	75	Black	36	Surgalloy™	VS-964
<b>KV-34</b> 1/2 Circle 37 mm 	1	4	2 × 75	White	12	D-Tach™	VS-766G

## SPATULA

Needles Specifications	Image	Size		Length cm	Color	Units Per box	Comments	Reorder code
		USP	EP					
<b>SE-100-8</b> 1/4 Circle 6.3 mm Double-arm Curvature 100°		8-0	0.4	45	Black	12	Surgalloy™	S-2176K
<b>SS-14</b> 1/4 Circle 8.6 mm Double-arm Curvature 115°		6-0	0.7	45	Black	12	Surgalloy™	S-1780K
<b>SS-24</b> 1/4 Circle 8.6 mm Double-arm Curvature 95°		6-0	0.7	45	Black	12	Surgalloy™	S-2780K
		4-0	1.5	45	White	12		S-2782K
<b>SE-140-8</b> 3/8 Circle 6.3 mm Double-arm Curvature 140°		8-0	0.4	30	Black	12	Surgalloy™	S-2752K
		8-0	0.4	45	Black	12	Surgalloy™	S-1734K
		7-0	0.5	45	Black	12	Surgalloy™	S-1733K
		6-0	0.7	45	Black	12	Surgalloy™	S-1732K

# V-LOC™ 90

Absorbable Wound Closure Device

Secure. Fast. Effective.†



**THE VLOC™ DEVICE**  
ADVANTAGE  
**NOW WITH STRONGER,  
SHARPER NEEDLES\***

Please refer to the individual package IFUs for complete instructions, indications, contraindications, warnings and precautions.



\* When compared with previous generation

† When compared with traditional suture



<b>Description</b>	The V-Loc™ 90 absorbable wound closure device consists of a barbed absorbable thread, armed with a surgical needle at one end and a loop end effector at the other. The barb and loop end effector design allow for tissue approximation without the need to tie surgical knots. The V-Loc™ 90 absorbable wound closure device is prepared from a synthetic polyester composed of glycolide, dioxanone, and trimethylene carbonate. The absorbable wound closure device is available undyed or violet. The device is sterile, inert, noncollagenous and nonantigenic.
<b>Indications</b>	V-Loc™ 90 absorbable wound closure devices are indicated for soft tissue approximation where use of an absorbable suture is appropriate. The safety and effectiveness of the V-Loc™ 90 absorbable wound closure device has not been established for use in fascial closures (abdominal wall, thoracic, extremity fascial closures), gastrointestinal anastomoses, cardiovascular anastomoses, neurological, ophthalmic, orthopedic or microsurgery.
<b>Actions</b>	V-Loc™ 90 absorbable wound closure devices elicit a minimal acute inflammatory reaction in tissues, which is followed by a gradual encapsulation of the device by fibrous connective tissue. Progressive loss of tensile strength and eventual absorption of V-Loc™ 90 absorbable wound closure device occurs by means of hydrolysis where the device is broken down to glycolic acid, dioxanoic acid, propane diol and carbon dioxide which are subsequently absorbed and metabolized by the body.
<b>Tensile Strength</b>	Approximately 90% of initial tensile strength at 1 week, 75% at 2 weeks postimplant.
<b>Sterilization</b>	Ethylene oxide.
<b>Absorbtion</b>	Absorption begins as a loss of tensile strength without appreciable loss of mass. Animal data indicates that the absorption is essentially complete between 90 and 110 days.
<b>Packaging</b>	The V-Loc™ 90 absorbable wound closure devices are available in USP and EP sizes 2-0 (3 metric) through 4-0 (1.5 Metric). The V-Loc™ 90 absorbable wound closure devices are available dyed violet or undyed (clear). The devices are supplied sterile, in pre-cut lengths affixed to various needle types. The devices are available in box quantities of one dozen. Please refer to the individual package IFUs for complete instructions, indications, contraindications, warnings and precautions.

## V-Loc™ 90 Absorbable

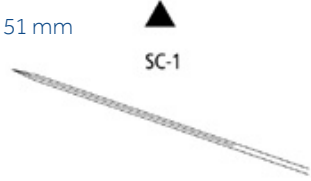
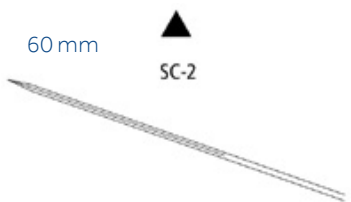
V-Loc™ 90 wound closure device is a uni-directional short-term absorbable knotless barbed suture.

PREMIUM REVERSE CUTTING								
Needles Specifications	Size USP	EP	Length cm	Color	Units Per box	Comments	Reorder code	
<b>P-11</b> 3/8 Circle 16 mm 	4-0	1.5	30	Undyed	12	Surgalloy™	VLOCM1413	
	4-0	1.5	45	Undyed	12	Surgalloy™	VLOCM1423	
<b>P-12</b> 3/8 Circle 19 mm 	4-0	1.5	15	Undyed	12	Surgalloy™	VLOCM0003	
	4-0	1.5	30	Undyed	12	Surgalloy™	VLOCM0013	
	4-0	1.5	45	Undyed	12	Surgalloy™	VLOCM0023	
	4-0	1.5	58	Undyed	12	Surgalloy™	VLOCM0033	
	3-0	2	15	Undyed	12		VLOCM0004	
	3-0	2	30	Undyed	12		VLOCM0014	
	3-0	2	45	Undyed	12		VLOCM0024	
	3-0	2	58	Undyed	12	Surgalloy™	VLOCM0034	
	2-0	3	15	Undyed	12		VLOCM0005	
	2-0	3	30	Undyed	12		VLOCM0015	
	2-0	3	45	Undyed	12		VLOCM0025	
	2-0	3	58	Undyed	12		VLOCM0035	




## PREMIUM REVERSE CUTTING

Needles Specifications	Size		Length cm	Color	Units Per box	Comments	Reorder code
	USP	EP					
<b>P-14</b> 3/8 Circle 24 mm 	<b>4-0</b>	1.5	30	Undyed	12	Surgalloy™	VLOCM0113
	<b>4-0</b>	1.5	45	Undyed	12	Surgalloy™	VLOCM0123
	<b>4-0</b>	1.5	58	Undyed	12	Surgalloy™	VLOCM0133
	<b>3-0</b>	2	30	Undyed	12	Surgalloy™	VLOCM0114
	<b>3-0</b>	2	45	Undyed	12	Surgalloy™	VLOCM0124
	<b>3-0</b>	2	58	Undyed	12	Surgalloy™	VLOCM0134
	<b>2-0</b>	3	30	Undyed	12		VLOCM0115
	<b>2-0</b>	3	45	Undyed	12		VLOCM0125
	<b>2-0</b>	3	58	Undyed	12		VLOCM0135
<b>P-17</b> 3/8 Circle 30 mm 	<b>3-0</b>	2	45	Undyed	12		VLOCM1824

## CONVENTIONAL CUTTING



Needles Specifications	Size		Length cm	Color	Units Per box	Comments	Reorder code
	USP	EP					
<b>SC-1</b> Straight 51 mm 	<b>3-0</b>	2	45	Undyed	12	Reliapoint™	VLOCM1624
<b>SC-2</b> Straight 60 mm 	<b>4-0</b>	1.5	45	Undyed	12	Reliapoint™	VLOCM0223
	<b>3-0</b>	2	45	Undyed	12		VLOCM0224

## TAPER POINT

Needles Specifications	Size USP	EP	Length cm	Color	Units Per box	Comments	Reorder code
<b>CV-15</b> 3/8 Circle 17 mm 	4-0	1.5	15	Violet	12		VLOCM1203
	3-0	2	15	Violet	12	Surgalloy™	VLOCM1204
<b>CV-23</b> 1/2 Circle 17 mm 	4-0	1.5	23		12		VLOCM0843
	4-0	1.5	30	Violet	12		VLOCM0813
	3-0	2	15	Undyed	12	Surgalloy™	VLOCM1904
	3-0	2	15	Violet	12	Surgalloy™	VLOCM0804
	3-0	2	23	Undyed	12	Surgalloy™	VLOCM1944
	3-0	2	23	Violet	12	Surgalloy™	VLOCM0844
	3-0	2	30	Violet	12		VLOCM0814
	3-0	2	45	Violet	12		VLOCM0824
<b>V-20</b> 1/2 Circle 26 mm 	4-0	1.5	15	Violet	12	Surgalloy™	VLOCM0603
	4-0	1.5	45	Violet	12	Surgalloy™	VLOCM0623
	3-0	2	15	Violet	12		VLOCM0604
	3-0	2	15	Undyed	12	Surgalloy™	VLOCM2004
	3-0	2	23	Violet	12	Surgalloy™	VLOCM0644
	3-0	2	23	Undyed	12	Surgalloy™	VLOCM2044
	3-0	2	30	Violet	12		VLOCM0614
	3-0	2	45	Violet	12		VLOCM0624
	2-0	3	45	Violet	12		VLOCM0625



## TAPER POINT

Needles Specifications	Size USP	EP	Length cm	Color	Units Per box	Comments	Reorder code
<b>GS-22</b> 1/2 Circle 27 mm 	2-0	3	15	Violet	12		VLOCM2105
	2-0	3	15	Undyed	12		VLOCM2205
	2-0	3	23	Violet	12		VLOCM2145
	2-0	3	23	Undyed	12		VLOCM2245
	2-0	3	30	Violet	12		VLOCM2115
<b>V-30</b> 1/2 Circle 30 mm 	3-0	2	23	Violet	12		VLOCM1544
	2-0	3	23	Violet	12		VLOCM1545
<b>GS-21</b> 1/2 Circle 37 mm 	3-0	2	23	Violet	12		VLOCM0344
	3-0	2	45	Violet	12		VLOCM0324
	2-0	3	15	Violet	12		VLOCM0305
	2-0	3	23	Violet	12		VLOCM0345
	2-0	3	30	Violet	12		VLOCM0315
	2-0	3	45	Undyed	12		VLOCM3225
<b>GU-46</b> 5/8 Circle 27 mm 	3-0	2	15	Violet	12	Surgalloy™	VLOCM1704
	3-0	2	23	Violet	12		VLOCM1744

# V-LOC™ 180

Absorbable Wound Closure Device

Secure. Fast. Effective.†



**THE VLOC™ DEVICE**  
ADVANTAGE  
**NOW WITH STRONGER,  
SHARPER NEEDLES\***

Please refer to the individual package IFUs for complete instructions, indications, contraindications, warnings and precautions.


\* When compared with previous generation



† When compared with traditional suture

<b>Description</b>	The V-Loc™ 180 absorbable wound closure device consists of a barbed absorbable thread, armed with a surgical needle at one end and a loop end effector at the other. The barb and loop end effector design allow for tissue approximation without the need to tie surgical knots. The V-Loc™ 180 absorbable wound closure device is prepared from a copolymer of glycolic acid and trimethylene carbonate. The absorbable wound closure device is available clear or green. The device is sterile, inert, noncollagenous and nonantigenic.
<b>Indications</b>	V-Loc™ 180 absorbable wound closure devices are indicated for soft tissue approximation where use of an absorbable suture is appropriate. The safety and effectiveness of the V-Loc™ 180 absorbable wound closure device has not been established for use in fascial closures (abdominal wall, thoracic, extremity fascial closures), gastrointestinal anastomoses, cardiovascular anastomoses, neurological, ophthalmic, orthopedic or microsurgery.
<b>Actions</b>	When implanted in animal tissue, V-Loc™ 180 absorbable wound closure devices elicit a minimal tissue reaction characteristic of foreign body response to substance. The tissue reaction resolves as the device is absorbed. Experimental animal studies for this absorbable material produced no evidence of antigenicity, systemic toxicity, carcinogenicity, mutagenicity, teratogenicity or adverse effects on reproductive performance.
<b>Tensile Strength</b>	Studies demonstrate the following approximate % of U.S.P. and E.P. knot pull minimum average specification for the labeled diameter size: Approximately 80% of initial tensile strength at 1 week, 75% at 2 weeks, and 65% at 3 weeks post implant.
<b>Sterilization</b>	Ethylene oxide.
<b>Absorbtion</b>	Animal data indicates that the absorption is minimal until about the 60th day after implantation. Absorption is substantially complete within six months.
<b>Packaging</b>	The V-Loc™ 180 absorbable wound closure devices are available in USP and EP sizes 4–0 (1.5 Metric) through 0 (3.5 Metric). The V-Loc™ 180 absorbable wound closure devices are dyed with D&C No. 6 or available Clear (Undyed). The devices are supplied sterile, in pre-cut lengths affixed to various needle types. The devices are available in box quantities of one dozen.


## V-Loc™ 180 Absorbable

V-Loc™ 180 wound closure device is a uni-directional long-term absorbable knotless barbed suture.

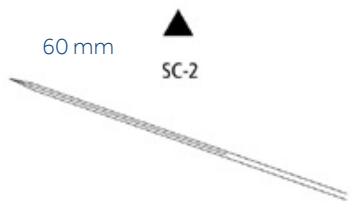
REVERSE CUTTING							
Needles Specifications	Size		Length cm	Color	Units Per box	Comments	Reorder code
	USP	EP					
<b>GS-11</b> 1/2 Circle    37 mm 	<b>2-0</b>	3	30	Clear	12	Reliapoint™	VLOCL2815
	<b>0</b>	3.5	45	Clear	12	Reliapoint™	VLOCL2826

PREMIUM REVERSE CUTTING							
Needles Specifications	Size		Length cm	Color	Units Per box	Comments	Reorder code
	USP	EP					
<b>P-11</b> 3/8 Circle    16 mm 	<b>4-0</b>	1.5	30	Green	12		VLOCL1413
<b>P-12</b> 3/8 Circle    19 mm 	<b>4-0</b>	1.5	15	Clear	12	Surgalloy™	VLOCL0003
	<b>4-0</b>	1.5	30	Clear	12	Surgalloy™	VLOCL0013
	<b>4-0</b>	1.5	45	Clear	12	Surgalloy™	VLOCL0023
	<b>4-0</b>	1.5	58	Clear	12	Surgalloy™	VLOCL0033
	<b>3-0</b>	2	15	Clear	12		VLOCL0004
	<b>3-0</b>	2	30	Clear	12		VLOCL0014
	<b>3-0</b>	2	45	Clear	12		VLOCL0024
	<b>3-0</b>	2	58	Clear	12		VLOCL0034
	<b>2-0</b>	3	30	Clear	12		VLOCL0015
	<b>2-0</b>	3	45	Clear	12		VLOCL0025


## PREMIUM REVERSE CUTTING

Needles Specifications	Size		Length cm	Color	Units Per box	Comments	Reorder code
	USP	EP					
<b>P-14</b> 3/8 Circle 24 mm 	4-0	1.5	45	Clear	12	Surgalloy™	VLOCL0123
	4-0	1.5	58	Clear	12	Surgalloy™	VLOCL0133
	3-0	2	30	Clear	12		VLOCL0114
	3-0	2	45	Clear	12		VLOCL0124
	3-0	2	60	Clear	12		VLOCL0134
	2-0	3	30	Clear	12		VLOCL0115
	2-0	3	45	Clear	12		VLOCL0125
	2-0	3	60	Clear	12		VLOCL0135






## CONVENTIONAL CUTTING

Needles Specifications	Size		Length cm	Color	Units Per box	Comments	Reorder code
	USP	EP					
<b>SC-2</b> Straight 60 mm 	4-0	1.5	45	Clear	12		VLOCL0223
	3-0	2	45	Clear	12		VLOCL0224


## TAPER POINT

Needles Specifications	Size		Length cm	Color	Units Per box	Comments	Reorder code
	USP	EP					
<b>CV-23</b> 1/2 Circle 17 mm 	4-0	1.5	15	Green	12		VLOCL0803
	4-0	1.5	30	Green	12	Surgalloy™	VLOCL0813
	3-0	2	15	Green	12		VLOCL0804
	3-0	2	15	Clear	12		VLOCL1904
	3-0	2	23	Green	12	Surgalloy™	VLOCL0844
	3-0	2	30	Green	12		VLOCL0814
	3-0	2	45	Green	12		VLOCL0824

## TAPER POINT

Needles Specifications	Size		Length cm	Color	Units Per box	Comments	Reorder code
	USP	EP					
<b>V-20</b> 1/2 Circle 26 mm 	4-0	1.5	15	Green	12		VLOCL0603
	4-0	1.5	30	Green	12		VLOCL0613
	3-0	2	15	Green	12		VLOCL0604
	3-0	2	23	Green	12	Surgalloy™	VLOCL0644
	3-0	2	30	Green	12	Surgalloy™	VLOCL0614
	3-0	2	45	Green	12	Surgalloy™	VLOCL0624
	2-0	3	30	Green	12		VLOCL0615
	2-0	3	45	Green	12		VLOCL0625
<b>GS-22</b> 1/2 Circle 27 mm 	2-0	3	15	Green	12		VLOCL2105
	2-0	3	23	Green	12		VLOCL2145
<b>V-30</b> 1/2 Circle 30 mm 	3-0	2	23	Green	12		VLOCL1544
	2-0	3	23	Green	12		VLOCL1545
<b>GS-21</b> 1/2 Circle 37 mm 	3-0	2	23	Green	12		VLOCL0344
	3-0	2	45	Green	12		VLOCL0324
	2-0	3	15	Green	12		VLOCL0305
	2-0	3	23	Green	12		VLOCL0345
	2-0	3	30	Green	12		VLOCL0315
	2-0	3	45	Green	12		VLOCL0325
	2-0	3	60	Green	12		VLOCL0335
	0	3.5	15	Green	12		VLOCL0306
	0	3.5	23	Green	12		VLOCL0346
	0	3.5	30	Green	12		VLOCL0316
	0	3.5	45	Green	12		VLOCL0326
	0	3.5	60	Green	12		VLOCL0336
<b>GS-25</b> 1/2 Circle 48 mm (Not to scale) 	0	3.5	30	Green	12		VLOCL0416
	0	3.5	45	Green	12		VLOCL0426
	0	3.5	60	Green	12		VLOCL0436

## TAPER POINT

Needles Specifications	Size		Length cm	Color	Units Per box	Comments	Reorder code
	USP	EP					
<b>SK</b> SKI                      23 mm 	<b>3-0</b>	2	23	Green	12		VLOCL2744

# V-LOC™ PBT

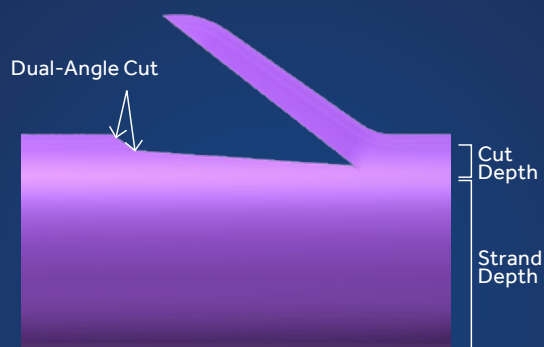
PBT Non Absorbable Wound Closure Device

Secure. Fast. Effective.†



## Dual-angle cut

Produces barbs with strong anchoring force while preserving the integrity of the strand's strength.



\* When compared with previous generation  
† When compared with traditional suture




Please refer to the individual package IFUs for complete instructions, indications, contraindications, warnings and precautions.





<b>Description</b>	The V-Loc™ Polybutester (PBT) non-absorbable wound closure device consists of a barbed non-absorbable thread, armed with a surgical needle at one end and a loop end effector at the other. The barb and loop end effector design allow for tissue approximation without the need to tie surgical knots. The V-Loc™ PBT nonabsorbable wound closure device is composed of a copolymer of butylene terephthalate and polytetramethylene ether glycol. The device is sterile, inert, noncollagenous and nonantigenic.
<b>Indications</b>	V-Loc™ PBT non-absorbable wound closure devices are indicated for soft tissue approximation. The safety and effectiveness of the V-Loc™ PBT non-absorbable wound closure device has not been established for use in fascial closures (abdominal wall, thoracic, extremity fascial closures), gastrointestinal anastomoses, cardiovascular anastomoses, neurological, ophthalmic, orthopedic or microsurgery.
<b>Actions</b>	Sofsilks™ sutures elicit a minimal acute inflammatory reaction in tissue, which is followed by gradual encapsulation of the suture by fibrous connective tissue.
<b>Tensile Strength</b>	No known significant change in tensile retention occurs in vivo for the V-Loc™ PBT non-absorbable wound closure device.
<b>Sterilization</b>	Gamma irradiation.
<b>Absorbtion</b>	Permanent.
<b>Packaging</b>	The V-Loc™ PBT non-absorbable wound closure devices are available in USP and EP sizes 1 (4 Metric) through 3–0 (2 Metric). The V-Loc™ PBT non-absorbable wound closure devices are pigmented with Copper Phthalocyanine Blue. The devices are supplied sterile, affixed to taperpoint needle types. The devices are available in box quantities of one dozen.

## V-Loc™ PBT

V-Loc™ PBT wound closure device is a uni-directional permanent knotless barbed suture.

TAPER POINT									
Needles Specifications			Size		Length cm	Color	Units Per box	Comments	Reorder code
			USP	EP					
<b>CV-23</b> 1/2 Circle 17 mm			<b>3-0</b>	2	30	Blue	12	Surgalloy™	VLOCN0814
<b>V-20</b> 1/2 Circle 26 mm			<b>3-0</b>	2	15	Blue	12		VLOCN0604
			<b>3-0</b>	2	23	Blue	12		VLOCN0644
			<b>3-0</b>	2	30	Blue	12	Surgalloy™	VLOCN0614
			<b>2-0</b>	3	15	Blue	12		VLOCN0605
			<b>2-0</b>	3	30	Blue	12		VLOCN0615
<b>GS-22</b> 1/2 Circle 27 mm			<b>0</b>	3.5	45	Blue	12		VLOCN1126

## TAPER POINT

Needles Specifications	Size		Length cm	Color	Units Per box	Comments	Reorder code
	USP	EP					
<b>GS-21</b> 1/2 Circle 37 mm 	<b>2-0</b>	3	15		12		VLOCN0305
	<b>2-0</b>	3	45	Blue	12		VLOCN0325
	<b>0</b>	3.5	15	Blue	12		VLOCN0306
	<b>0</b>	3.5	45	Blue	12		VLOCN0326
	<b>1</b>	4	45	Blue	12		VLOCN0327
<b>GS-25</b> 1/2 Circle 48 mm 	<b>1</b>	4	60	Blue	12		VLOCN0437

# PLEDGETS

Non-Absorbable, undyed fabrics constructed of PTFE

**Simple strength. Secure placement.**



**Compared to rectangular pledgets:**



Oval pledgets conform better to the shape of the anatomical structure and result in less foreign mass in the body<sup>2</sup>



The lack of corners prevents the oval pledget from catching the wrong orientation and the edges are less likely to overlap<sup>2</sup>



Soft oval pledgets are compliant, yet not fuzzy<sup>2</sup>

**Description**

PTFE Pledgets are non-absorbable, undyed fabrics composed of poly (tetrafluoroethylene).

**Indications**

PTFE Pledgets are indicated for use as non-absorbable suture buttresses when there is a possibility of sutures tearing through tissue. These pledgets are used in various surgical procedures such as vascular closure, septal repair, myocardial closure, valvular suturing and hepatic repair.

**Actions**

PTFE Pledgets elicit a minimal acute inflammatory reaction in tissues, followed by gradual encapsulation of the fabric by fibrous connective tissue.

**Absorption**

PTFE Pledgets are not absorbed, nor is any significant change in tensile strength retention known to occur in vivo.

**Packaging**

PTFE Pledgets are available sterile, prepunched, square, rectangular or oval in shape and in box quantities of 3 dozen.

# PLEDGETS

Our broad range of pledgets are offered as free (bulk) and on sutures.

## How to Read a Suture Code

**Pledgets with suture example: VCD20316MFE**

**V**= Vascular, **CD**= Coated Dacron, **20**= Size 2-0, **316**= Needle, **M**= Multipack, **F**= Firm (**S**=Soft), **E**= Ellipse

**Free pledget example: VPG36FE**

**V**= Vascular, **PG**= Pledgets, **36**= 3x6 mm (**56**= 4.5x6 mm), **F**= Firm (**S**=Soft), **E**= Ellipse

OVAL PLEDGETS	
Description	Order Code
Ti-Cron™ 2-0 BL/WH 8x30 CV-316 DA 3x6F	VCD20316MFE
Ti-Cron™ 2-0 BL/WH 8x30 CV-316 DA 3x6S	VCD20316MSE
Ti-Cron™ 2-0 BL 30 CV-316 DA 3x6F	VCD20316FE
Ti-Cron™ 2-0 BL 30 CV-316 DA 3x6S	VCD20316SE
Ti-Cron™ 2-0 BL/WH 8x30 CV-305 DA 3x6F	VCD20305MFE
Ti-Cron™ 2-0 BL/WH 8x30 CV-305 DA 3x6S	VCD20305MSE
Ti-Cron™ 2-0 BL 30 CV-305 DA 3x6F	VCD20305FE
Ti-Cron™ 2-0 BL 30 CV-305 DA 3x6S	VCD20305SE
Ti-Cron™ 2-0 BL/WH 8x30 Y31 DA 3x6F	VCD2031MFE
Ti-Cron™ 2-0 BL/WH 8x30 Y31 DA 3x6S	VCD2031MSE
Pledget Oval 3x6 Firm 10 PK	VPG36FE
Pledget Oval 3x6 Soft 10 PK	VPG36SE
Pledget Oval 4.5x6 Firm 10 PK	VPG56FE
Pledget Oval 4.5x6 Soft 10 PK	VPG56SE
Ti-Cron™ 2-0 BL/WH 8x30 CV-331 DA 3x6F	VCD20331FE
Ti-Cron™ 2-0 BL/WH 8x36 CV-305 DA 3x6F	VCD2036305FE
Ti-Cron™ 2-0 BL/WH 8x30 CV-305 DA 4.5F	VCD2030556MFE
Ti-Cron™ 2-0 BL/WH 8x30 CV-305 DA 4.5S	VCD2030556MSE

RECTANGURAL & SQUARE PLEDGETS	
Description	Order Code
Ti-Cron™ 2-0 5 Blue 5 White 30" CV-316 DA 3x7F	8886332356
Ti-Cron™ 2-0 5 Blue 5 White 30" CV-305 DA 3x7F	8886332456
Ti-Cron™ 2-0 4 Blue 4 White 30" Y-31 DA 3x7F	8886321956
Ti-Cron™ 2-0 3 Blue 3 White 36" CV-305 DA 3x7F	8886330656
Ti-Cron™ 2-0 4 Blue 4 White 30" CV-331 DA 3x7F	8886309156
Ti-Cron™ 2-0 4 Blue 4 White 30" Y-5 DA 3x7F	8886321356
Ti-Cron™ 2-0 4 Blue 4 White 30" Y-31 DA 3x7S	VCD656
Ti-Cron™ 2-0 4 Blue 4 White 30" CV-331 DA 3x3F	8886308656
Ti-Cron™ 2-0 4 Blue 4 White 30" Y-31 DA 3x3F	8886321856
Ti-Cron™ 2-0 5 Blue 5 White 30" CV-305 DA 4.8x9.5F	8886303756
Ti-Cron™ 2-0 4 Blue 4 White 30" Y-31 DA 3x5F	VCD110
3x7 Firm Pledget	8886867701

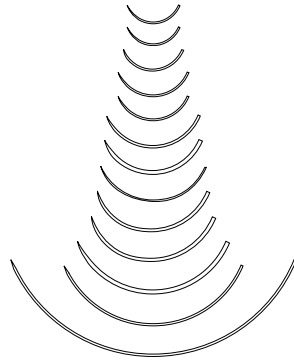
# WOUND CLOSURE ACCESSORIES

WOUND CLOSURE ACCESSORIES	
Description	Reorder code
Umbilical Tape	88868619-03
Polyester Tape	88868618-00
Bone-wax	BW25G
Bolsters	88868631-12

# NEEDLE CHART

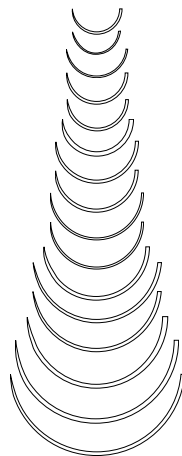
## CARDIOVASCULAR

Code	Comp	Length
⊙ taper point - $\frac{3}{8}$ circle		
CV-1/CVF-1	BV-1	9 mm
CV-310	BV	10 mm
CV/CVF	BV	11 mm
CV-11/CVF-11	C-1	13 mm
CV-301	C-1*	13 mm
CV-307	BB*	16 mm
CV-15/CVF-15	BB	17 mm
CV-327	BB*	18 mm
CV-16		20 mm
CV-13	BB-1	22 mm
CV-18	TE-2	26 mm
CV-17	TE-2	32 mm
CV-19		50 mm



### ⊙ taper point - $\frac{1}{2}$ circle

CV-20	RB-3	10 mm
CV-337	RB-3	10 mm
CVF-21	TF	12 mm
CV-22/CVF-22	RB-2	13 mm
CV-330	RB-2	13 mm
CV-26	RB-1*	15 mm
CV-331	RB-1*	16 mm
CV-23/CVF-23	RB-1	17 mm
CV-316	SH-2	20 mm
CV-24	SH-2	20 mm
CV-25	SH-1	22 mm
CV-305	SH*	25 mm
V-20/VF-20	SH	26 mm
V-30		30 mm
CV-300	MH*	35 mm
V-26	MH*	37 mm



## CARDIOVASCULAR

Code	Comp	Length
⊙ taper point - CardioCurve™		
CVL		11 mm
CVL-11		13 mm



### ⊙ penetrating taper - $\frac{3}{8}$ circle

YE-1	CC-1	13 mm
YE-7	V-4	18 mm



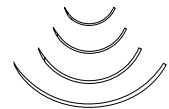
### ⊙ penetrating taper - $\frac{1}{2}$ circle

Y-31	V-5	16 mm
Y-16	V-6*	20 mm
Y-5	V-7	26 mm



### ⊙ tapercutting - $\frac{3}{8}$ circle

KV-1/KVF-1	CC	9 mm
KV-11/KVF-11	CC-1	13 mm
KV-15	V-4	17 mm
KV-26	V-26	26 mm



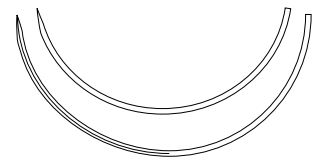
### ⊙ tapercutting - $\frac{1}{2}$ circle

KV-5/KVF-5	V-5	17 mm
KV-7	V-7	26 mm



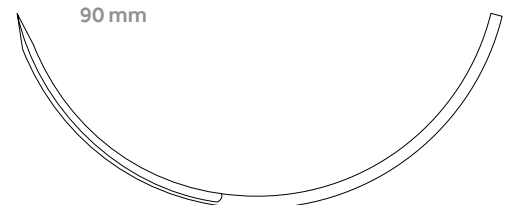
### ▲ conventional cutting - $\frac{1}{2}$ circle

SCC	CCS	48 mm
SCC-1	CCS-1	60 mm



### ▲ conventional cutting - $\frac{3}{8}$ circle

GCC-90	90 mm
--------	-------





**PLASTIC/COSMETIC**

Code	Comp	Length
------	------	--------

▲ premium conventional cutting - 3/8 circle

PC-13	P-1	11 mm	
PC-10	PC-1	13 mm	
PC-11	PC-3	16 mm	
PC-12	PC-5	19 mm	

▼ premium reverse cutting - 3/8 circle

P-16	P-6	7 mm	
P-10	P-1	11 mm	
P-13	P-3	13 mm	
P-15		15 mm	
P-11	PS-3	16 mm	
P-12	PS-2	19 mm	
P-14	PS-1	24 mm	
P-17	PSL	30 mm	
P-18	PSLX	39 mm	

▼ premium reverse cutting - 1/2 circle

P-21	P-2	9 mm	
P-22	PS-5	13 mm	
P-24	PS-4*	16 mm	

◊ DermaX™ premium X cutting · 3/8 circle

DX-11		11 mm	
DX-13		13 mm	
DX-16		16 mm	
DX-19		19 mm	

◊ DermaX™ premium X cutting · 1/2 circle

DXH-16		16 mm	1/2
--------	--	-------	-----

**CUTICULAR**

Code	Comp	Length
------	------	--------

▼ reverse cutting - 3/8 circle

C-1	C-2*	12 mm	
C-12	FS-3	16 mm	
C-13	FS-2	19 mm	
C-14	FS-1	24 mm	
C-15	FS	26 mm	
C-16	FSL	30 mm	
C-18		35 mm	
C-17	FSLX*	39 mm	
C-50	LS-1*	50 mm	
GS-18	LS-1*	77 mm	

▼ reverse cutting - 1/2 circle

C-21	M-1*	13 mm	
C-26	J-1*	15 mm	
C-22	J-1	18 mm	
C-27		20 mm	
C-23	X-1*	24 mm	
C-25	CP-2*	30 mm	

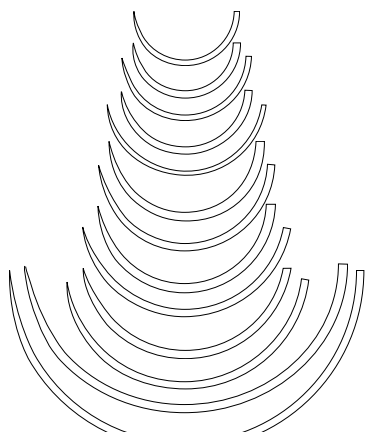
▲ cutting - straight

SC-11	SC-1*	11 mm	
SC-4		24 mm	
SC		44 mm	
SC-1	TS*	51 mm	
SC-2	KS	60 mm	
SC-250		254 mm	





\*Suggested

## GENERAL SURGERY

Code	Comp	Length
⊙ taper point - ½ circle		
GS-23	CT-3	22 mm
HGS-23	MO-7	22 mm
GS-22	CT-2	27 mm
HGS-22	MO-6	27 mm
GS-20	CT-1*	34 mm
HGS-20	MO-5*	34 mm
GS-21	CT-1	37 mm
HGS-21	MO-4	37 mm
GS-24	CT	40 mm
HGS-24	MO-2	40 mm
GS-25	CTX	48 mm
GS-26	TP-1	65 mm
GS-27	XLH	76 mm



### ⊙ taper point - straight

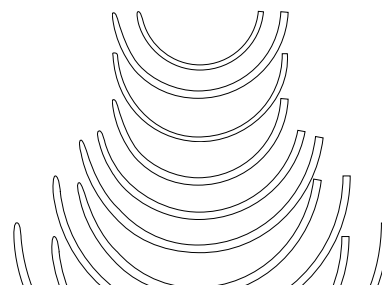
ST-4	ST-4	19 mm	
ST-3		38 mm	
ST-2		44 mm	
ST-1	ST	51 mm	
ST		60 mm	

### ⊙ taper point - ski

SK		23 mm	
----	--	-------	--

### ● blunt taper point - Protect•Point™ - ½ circle

BTV-20	SHB	26 mm
BGS-21	CTB-1	37 mm
HBSG-21	MOB-4	37 mm
BTP-1	CTB-1	37 mm
BGS-24		40 mm
BGS-25	CTX-B	48 mm
BTP-X	CTX-B	48 mm
BGS-29	BP-1	65 mm
BGST-29	BP-1	65 mm
BGS-28	BP*	85 mm



### ● blunt point - Protect•Point™ - ½ circle

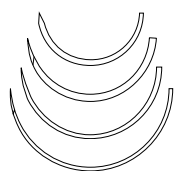
BP-27		64 mm	
-------	--	-------	---

### ⊖ diamond point - ⅓ circle

KV-20		20 mm	
-------	--	-------	---

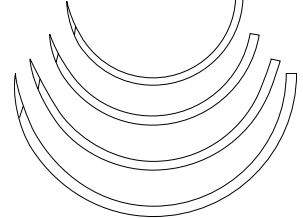
## GENERAL SURGERY

Code	Comp	Length
⊖ diamond point - ½ circle		
KV-25		22 mm
KV-30		27 mm
KV-8		30 mm
KV-9		35 mm




### ⊖ tapercutting - ½ circle

KV-34	V-34	37 mm
KV-37	V-37	40 mm
KV-40	V-40	48 mm
KV-56	V-56	60 mm



## MICROSURGERY

Code	Comp	Length
⊙ taper point - ⅓ circle		
MV-100-3	BV-100-3	4 mm
MV-50-3	BV-50-3	4 mm
MV-70-3	BV-75-3	4 mm
MV-135-3	BV-130-3*	4 mm
MV-70-4	BV-75-4	5 mm
MV-135-4	BV-130-4*	5 mm
MV-100-4	BV-100-4	5 mm
CV-345	BV-130-5	6 mm
MVF-135-5	BV-130-5	6 mm
MV-135-5	BV-130-5	6 mm
CV-351	BV-175-6	8 mm
MV-175-8	BV-175-6	8 mm
MVF-175-8	BV-175-6	8 mm
MV-175-9	BV-175-8	10 mm



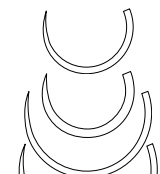
### ⊖ tapercutting - ⅓ circle

MVK-70-3	BV-75-3	4 mm
MVK-100-4	BV-100-4	5 mm



## UROLOGY

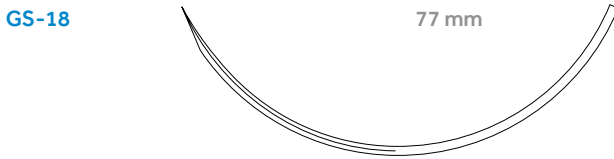
Code	Comp	Length
⊙ taper point - ⅓ circle		
GU-46	UR-6	27 mm
HGU-46	UR-6	27 mm
GU-45	UR-5	37 mm
GU-44	UR-4	40 mm



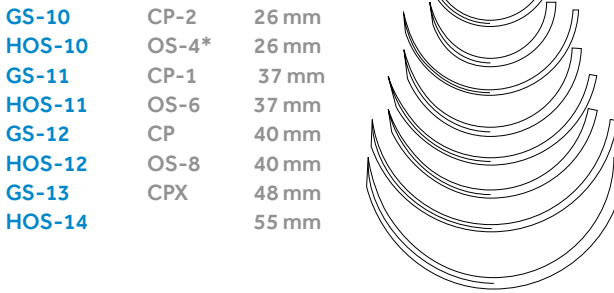
## ORTHOPEDIC

Code	Comp	Length
------	------	--------

### ▼ reverse cutting - 3/8 circle



### ▼ reverse cutting - 1/2 circle



### 🔪 tapercutting - straight



\*Suggested

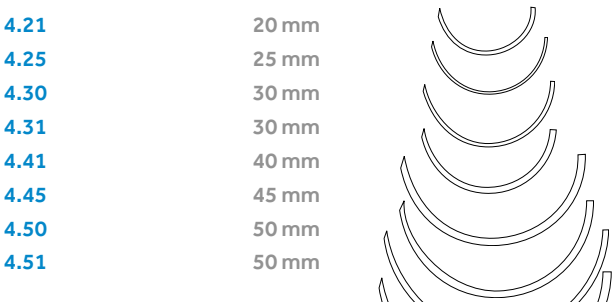
## MULTI-PURPOSE

Code	Comp	Length
------	------	--------

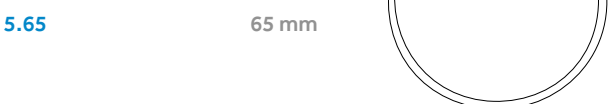
### 🔪 multi-purpose - 3/8 circle



### 🔪 multi-purpose - 1/2 circle



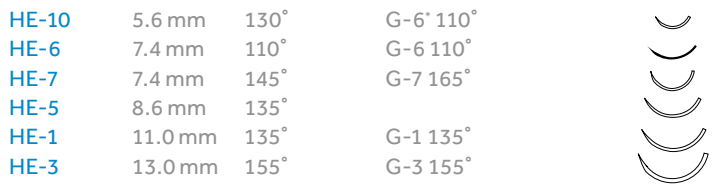
### 🔪 multi-purpose - 3/4 circle



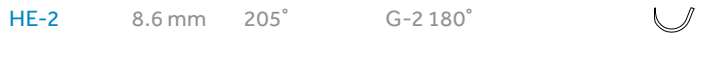
## OPHTHALMIC

Code	Length	Curvature	Comp
------	--------	-----------	------

### ▼ reverse cutting - 3/8 circle



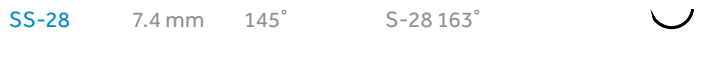
### ▼ reverse cutting - 1/2 circle



### 🔪 inverted premium spatula - 1/4 circle



### 🔪 inverted premium spatula - 3/8 circle



### 🔪 premium spatula - 1/8 circle



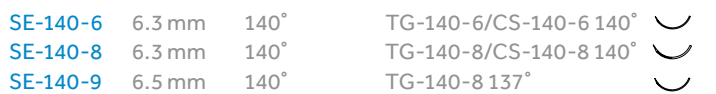
### 🔪 premium spatula - 1/4 circle



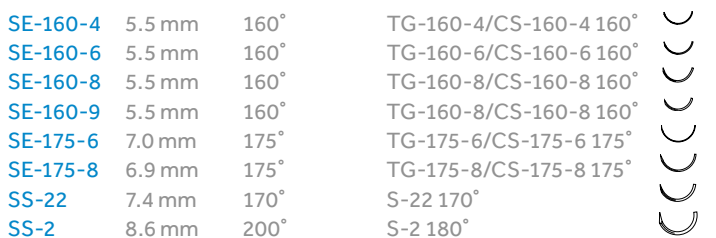
### 🔪 premium spatula - 5/16 circle



### 🔪 premium spatula - 3/8 circle



### 🔪 premium spatula - 1/2 circle



### 🔪 premium spatula - compound curve



# ORDERING INFORMATION

<b>Velosorb™</b>		CV-485	<b>10</b>	SV-9944G	<b>13</b>
CV-8977	<b>8</b>	CV-486	<b>10</b>	CV-9916G	<b>13</b>
CV-8643	<b>8</b>	SV-9913	<b>11</b>	CV-1213	<b>13</b>
CV-8978	<b>8</b>	SV-32	<b>11</b>	CV-1214	<b>13</b>
CV-8644	<b>8</b>	SV-4893	<b>11</b>	CV-1215	<b>13</b>
CV-8979	<b>8</b>	SV-9915	<b>11</b>	CV-3100	<b>13</b>
CV-8653	<b>8</b>	SV-493	<b>11</b>	CV-3110	<b>13</b>
CV-8980	<b>8</b>	SV-494	<b>11</b>	CV-3120	<b>13</b>
SV-2289	<b>8</b>	SVP-494	<b>11</b>	CV-415	<b>14</b>
SV-1290	<b>8</b>	SVP-2296	<b>11</b>	CV-416G	<b>14</b>
SV-2290	<b>8</b>	SVP-2294	<b>11</b>	CV-417G	<b>14</b>
SV-2296	<b>8</b>	SVP-2295	<b>11</b>	CV-418	<b>14</b>
SV-2294	<b>8</b>	SV-496G	<b>11</b>	CV-55	<b>14</b>
SV-2295	<b>8</b>	SV-496	<b>11</b>	CV-56	<b>14</b>
SV-2910	<b>9</b>	SV-426	<b>11</b>	CV-57	<b>14</b>
SV-2920	<b>9</b>	SV-412	<b>11</b>	CV-318	<b>14</b>
SV-2297	<b>9</b>	SV-497G	<b>11</b>	CV-319	<b>14</b>
SV-2930	<b>9</b>	SV-425	<b>11</b>	CV-9975G	<b>14</b>
SV-2298	<b>9</b>	SV-424	<b>11</b>	CV-322	<b>14</b>
SV-2299	<b>9</b>	SVP-2251	<b>11</b>	CV-944	<b>14</b>
SV-2251	<b>9</b>	SV-6830	<b>11</b>	CV-323	<b>14</b>
SV-4431	<b>9</b>	SV-935	<b>11</b>	CV-945	<b>14</b>
SV-2252	<b>9</b>	SV-6840	<b>11</b>	CV-1324	<b>14</b>
SV-2253	<b>9</b>	SV-439	<b>11</b>	CV-946	<b>14</b>
SV-9940G	<b>9</b>	SVP-2285	<b>11</b>	CV-947	<b>14</b>
SV-9941G	<b>9</b>	SV-9910G	<b>12</b>	CV-9592	<b>15</b>
SV-2285	<b>9</b>	SV-9911G	<b>12</b>	CV-9591	<b>15</b>
SV-2286	<b>9</b>	SV-834G	<b>12</b>	CV-9590	<b>15</b>
SV-2287	<b>9</b>	SV-835G	<b>12</b>	CV-606	<b>15</b>
SV-9961G	<b>9</b>	SV-489G	<b>12</b>	CV-604	<b>15</b>
SV-4550	<b>10</b>	SV-844G	<b>12</b>	CV-605	<b>15</b>
SV-9926	<b>10</b>	SV-845G	<b>12</b>	UV-412	<b>15</b>
SV-9927	<b>10</b>	SV-488G	<b>12</b>	UV-401	<b>15</b>
SV-9928	<b>10</b>	SV-487G	<b>12</b>	UV-1381	<b>15</b>
CV-9936G	<b>10</b>	SV-9909G	<b>12</b>	UV-1382	<b>15</b>
CV-931G	<b>10</b>	SV-2238	<b>13</b>	CV-1673	<b>15</b>
CV-4721	<b>10</b>	SV-2239	<b>13</b>	CV-2293	<b>16</b>
CV-932G	<b>10</b>	SV-2240	<b>13</b>	CV-2275G	<b>16</b>
CV-4731	<b>10</b>	SV-9943	<b>13</b>	CV-2279	<b>16</b>

CV-2280	16	L-14	20	L-1551-K	23
CV-1674G	16	L-34	20	SL-171	23
CV-2003G	16	L-404	20	SL-172	23
CV-1002	16	L-40	20	L-4323K	23
CV-916	16	L-75	20	SL-630	24
CV-917	16	L-5	20	SL-607	24
CV-927	16	L-35	20	SL-605	24
CV-919	16	L-25	20	SL-608	24
CV-918	16	L-15	20	SL-639	24
CV-0146	16	L-36	20	SL-634	24
CV-2218	17	L-41	21	SL-690	24
CV-2217	17	L-405	21	SL-638	24
CV-2216	17	L-76	21	SL-691G	24
CV-0046	17	L-28	21	SL-691	24
CV-7192	17	L-26	21	SL-635G	24
CV-7191	17	L-42	21	SL-635	24
CV-7190	17	L-406	21	SL-637	24
CV-9953	17	L-111	21	SL-822	24
CV-501G	17	L-102	21	SL-822G	24
CV-500G	17	L-112	21	SL-636	24
		L-103	21	SL-695	24
<b>Polysorb™</b>		L-113	21	SL-692	24
L-71	20	L-114	21	SL-693	24
L-20	20	L-104	21	SL-694	24
L-11	20	L-116	21	SL-823	24
L-1-	20	L-115	21	CL-823	24
L-72	20	LL-111	22	SL-73-MG	24
L-22	20	LL-101	22	SL-586	24
L-2-	20	LL-221	22	SL-587	24
L-12	20	LL-112	22	SL-696	25
L-32	20	LL-102	22	SL-714	25
L-402	20	LL-322	22	SL-611	25
L-38	20	LL-222	22	SL-713	25
L-73	20	LL-113	22	SL-610	25
L-3	20	LL-103	22	SL-712	25
L-23	20	LL-323	22	CL-976	25
L-13-	20	LL-223	22	CL-977	25
L-33	20	LL-233	22	SL-674	25
L-43	20	LL-114	22	SL-675	25
L-403	20	LL-104	22	CL-85	25
L-39	20	LL-224	22	SL-85-MG	25
L-74	20	LL-115	22	SL-83-MG	25
L-4	20	LL-225	22	SL-632	25
L-27	20	LL-116	22	SL-613	25
L-24	20	LL-226	22	SL-631	25

SL-612	25	CL-528	27	SL-5636	29
CL-53MG	25	CL-833	27	SL-5629	29
CL-868	25	CL-595	27	SL-5627-G	29
CL-52-MG	25	CL-804	27	SL-5627	29
CL-869	25	CL-596	27	SL-5637	29
CL-769	25	CL-933	27	SL-5628G	29
CL-870	25	CL-805	27	SL-5628	29
CL-867	25	CL-597	27	SL-5638	29
CL-356-MG	26	CL-936	27	SL-5641G	30
CL-354-MG	26	CL-937	27	SL-5641	30
CL-55-MG	26	CL-382-MG	27	SL-5678	30
CL-355-MG	26	CL-380-MG	27	SL-5640G	30
CL-358-MG	26	CL-376-MG	27	SL-5640	30
CL-581	26	CL-381-MG	27	SL-5679G	30
CL-547	26	CL-377-MG	27	SL-5679	30
CL-552	26	CL-97-MG	27	SL-5680	30
SL-622	26	CL-383-MG	27	SL-1596	30
SL-721	26	CL-568	28	SL-1510G	30
SL-722	26	CL-554	28	SL-1644G	30
CL-815	26	CL-927	28	SL-1645G	30
CL-542	26	CL-569	28	SL-1613	30
CL-816	26	CL-555	28	SL-1614	30
CL-543	26	CL-928	28	SL-1614G	30
CL-916	26	CL-570	28	SL-1615	30
CL-817	26	CL-939	28	SL-5632G	30
CL-544	26	CL-929	28	SL-5632	30
CL-917	26	CL-979	28	SL-5633G	30
CL-534	26	SL-1624G	28	SL-5633	30
CL-818	26	SL-1625	28	SL-1660	31
CL-545	26	SL1625G	28	SL-1659	31
CL-918	26	SL-1689	28	SL-1658	31
CL-535	26	SL-1689G	28	SL-1655G	31
CL-919	26	SL-5688	29	SL-1656	31
CL-533	26	SL-5688-G	29	SL-1657	31
CL-982	26	SL-5687G	29	SL-1653	31
CL-992	26	SL-5687	29	SL-1654	31
CL-374-MG	27	SL-5691G	29	SL-654	31
CL-369-MG	27	SL-5690	29	SL-644	31
CL-526	27	SL-5626	29	SL-653	31
CL-536	27	SL-5686	29	SL-643	31
CL-375-MG	27	SL-5686G	29	SL-652	31
CL-370-MG	27	SL-5620	29	SL-642	31
CL-537	27	SL-5621	29	L-2800	32
CL-538	27	SL-5636-G	29	GL-889-2	32

GL-890	32	EL-440	34	CL-831	36
GL-881	32	GL-32-MG	34	CL-330	36
GL-882	32	GL-183	34	CL-885	36
CL-304	32	GL-282	34	CL-885-M	36
GL-303	32	CL-892	34	CL-832	36
GL-302	32	CL-42-MG	34	CL-331	36
CL-457	32	CL-893	34	CL-232	36
GL-522	32	CL-894	34	GL-222	36
GL-523	32	GL-120	35	GL-223	36
GL-524	32	GL-64-MG	35	GL-228	36
UL-102	33	GLS-321	35	GL-223-M	36
UL-101	33	GLS-121	35	GL-224	36
SL-432	33	GL-68-M	35	GL-225	36
SL-431	33	GL-321	35	CL-809	37
GL-45-MG	33	GL-121	35	CL-8-MG	37
GL-885	33	GL-69-MG	35	CL-842	37
GL-44-MG	33	GL-63-MG	35	CL-810	37
GL-884	33	GLS-322	35	CL-922	37
CL-782	33	GLS-122	35	CL-952	37
UL-201-2	33	GL-67-M	35	CL-9-MG	37
UL-202	33	GL-322	35	CL-13-MG	37
UL-213	33	GL-122	35	CL-843	37
GL-46-MG	33	GL-226	35	CL-811	37
UL-303	33	GL-62-MG	35	CL-4M	37
UL-212	33	GL-61-MG	35	CL-923	37
UL-47-MG	33	GLS-323	35	CL-953	37
GL-47-MG	33	GLS-123	35	CL-14-MG	37
UL-204	33	GL-66-M	35	CL-10-MG	37
UL-215	33	GL-323	35	CL-844	37
UL-216	33	GL-123	35	CL-812	37
UL-205	33	GL-125	35	CL-5-M	37
UL-214	33	GL-324	35	CL-2M	37
UL-203	33	GL-124	35	CL-812	37
CL-219	33	GL-332	35	CL-954	37
CL-460	34	GL-333	35	CL-924	37
CL-461	34	CL-882	36	CL-15-MG	37
CL-462	34	CL-829	36	CL-11-MG	37
GL-34-MG	34	CL-16-MG	36	CL-845	37
GL-181	34	CL-830	36	CL-813	37
GL-191	34	CL-883	36	CL-6-M	37
GL-33-MG	34	CL-17-MG	36	CL-925	37
GL-182	34	CLS-884	36	CL-955	37
GL-192	34	CL-886-M	36	CL-814	37
GL-182-M	34	CL-884	36	CL-926	37

GL-126	38	CL-665-M	39	CL-677	42
GL-127	38	CL-906	39	UL-676	42
GL-128	38	CL-63	39	CL-510	42
GL-129	38	CL-60-M	39	CL-511	42
CL-800	38	CL-968	39	CL-942	42
CL-962	38	CL-824	39	CL-837	42
CL-912	38	CL-659-M	39	CL-943	42
CL-22-MG	38	CL-969	39	CL-946	42
CL-853	38	CL-880	39	CL-838	42
CL-801	38	CL-825	39	CL-515	42
CL-913	38	CL824L	40	CL-944	42
CL-963	38	CL-835-L	40	CL-839	42
CL-20-MG	38	CL-53	40	CL-516	42
CL-802	38	SL-603	40	CL-945	42
CL-914	38	UL-879	40	CL-517	42
CL-964	38	UL-878	40	CL-603	43
CL-21-MG	38	UL-877	40	CL-602	43
CL-803	38	UL-876	40	CL-739	43
CL-855	38	UL-240	40	CL-840	43
CL-915	38	UL-245	40	CL-940	43
CL-965	38	UL-246	40	CL-841	43
CL-623-M	38	UL-146	40	CL-941	43
CL-622-M	38	EL415	41	CL-740-L	43
CL-779	38	EL450	41	CL-741-L	43
CL-911	38	CL-673G	41	GLT-69-MG	43
CL-990	38	CL-673	41	CL-956	43
CL-862	39	CL-86-MG	41	CL-259	43
CL-972	39	CL-87-MG	41	CL-250-M	43
CL-36-MG	39	CL-871	41	CL-240-M	43
CL-863	39	CL-996-MG	41	CL-957	43
CL-973	39	CL-997-MG	41	CL-958	43
CL-903	39	CL-998-MG	41	CLT-952	44
CL-35-MG	39	CL-202-M	41	CLT-9-MG	44
CL-30-MG	39	CL-930	41	CLT-13-MG	44
CL-864	39	CL-970	41	CLT-259	44
CL-974	39	CL-200-M	41	CLT-953	44
CL-904	39	CL-931	41	CLT-250-MG	44
CL-31-MG	39	CL-71-M	41	CLT-957	44
CL-865	39	CL-971	41	CLT-954	44
CL-64-M	39	CL-201-M	41	CLT-15-MG	44
CL-975	39	CL-932	41	CLT-11-MG	44
CL-905	39	CL-239	41	CLT-958	44
CL-633-M	39	CL-934	41	CLT-955	44
CL-644-M	39	CL-871-L	42	CL-950	44
CL-866	39	CL-576	42	CL-951	44



CLT-850-MG	44	L-2749K	49	SC-643	55
CLT-20-MG	44	L-1787K	49	XC-834	55
CLT-851-MG	44	L-0223K	49	XC-844	55
CLT-21-MG	44			XC-845	55
CLT-855	44	<b>Caprosyn™</b>		XC-823	56
CLT-623-M	44	C-212	52	UC-202	56
CL-948	45	C-214	52	UC-203	56
CL-949	45	SC-3634	52	UC-214	56
2-CL-852	45	SC-690G	52	UC-213	56
CL-271	45	SC-691	52	UC-204	56
CL-270	45	SC-694G	52	UC-403	56
CLT-633-M	45	SC-692	52	GC-33-M	56
CL-513	46	SC-693	52	UC-404	56
CL-730	46	SC-693G	52	UC-863	56
CL-005	46	SC-587	52	GC-321	57
CL-015	46	SC-586	52	GC-121	57
CL-030	46	SC-585	52	GC-122	57
CL-045	47	CC-543	53	GC-322	57
UL-046	47	CC-983	53	GC-123	57
CL-047	47	SC-3625G	53	GC-323	57
CL-055	47	SC-3626G	53	GC-124	57
CL-060	47	SC-5688G	53	CC-883	57
L-1750K	47	SC-5588G	53	CC-830	57
L-1756K	47	SC-5689G	53	CC-884	57
L-1755K	47	SC-5587G	53	CC-831	57
L-1791K	48	SC-5690G	53	CC-885	57
L-1670-K	48	SC-5590G	53	CC-810	57
L-1792K	48	SC-5635G	54	CC-922	57
L-1669-K	48	SC-5626G	54	CC-843	57
L-1740-K	48	SC-5626	54	CC-811	57
L-1742K	48	SC-5686G	54	CC-923	57
L-1745K	48	SC-5616G	54	CC-844	57
L-2752K	48	SC-5617G	54	CC-812	57
L-1739-K	48	SC-5627G	54	CC-924	57
L-1785K	48	SC-5637-G	54	CC-845	57
L-2753K	48	SC-5618G	54	CC-813	57
L-2750K	48	SC-5638-G	54	CC-925	57
L-2751K	48	SC-5679	54	CC-801	58
L-2747K	48	SC-5779	54	CC-915	58
L-2748K	48	SC-5640	54	CC-802	58
L-5118K	48	SC-5780G	54	CC-913	58
L-1796K	48	SC-5643G	54	CC-803	58
L-1897K	48	SC-543	55	CC-914	58
L-1562-K	48	SC-542	55	CC-872	58
L-2574K	49	SC-644	55	CC-904	58

CC-975	58	SM-5637	64	GM-323	68
CC-905	58	SM-5628	64	GM-125	68
UC-879	58	SM-5638	64	GM-421	68
UC-978	58	SM-5679	65	GM-432	68
CC-658	58	SM-5678	65	GM-332	68
CC-1945G	59	SM-5679G	65	CM-882	68
CC-196	59	CM-401	65	CM-883	68
		CM-400	65	CM-830	68
<b>Biosyn™</b>		SM-644	65	CM-884	68
M-31-	62	SM-643	65	CM-831	68
M-72	62	XM-834	66	CM-832	68
M-32-	62	XM-845	66	GM-222	68
M-33-	62	XM-822	66	GM-223	68
M-23-	62	XM-823	66	GM-326	68
M-74	62	SM-432	66	GM-327	68
M-75	62	SM-431	66	GM-226	69
M-114	62	SM-434	66	GM-227	69
SM-690	62	CM-126	66	CM-738-M	69
SM-791	62	UM-986	67	CM-842	69
SM-691G	62	UM-886	67	CM-810	69
SM-691	62	GM-875	67	CM-922	69
SM-822	62	UM-201	67	CM-9M	69
SM-922	62	UM-202	67	CM-843	69
SM-693	63	UM-213	67	CM-811	69
SM-923	63	UM-214	67	CM-953	69
SM-823	63	UM-203	67	CM-923	69
SM-633	63	UM-215	67	CM-14-M	69
SM-83-M	63	UM-204	67	CM-10-M	69
SM-632	63	UM-403	67	CM-844	69
CM-879	63	CM-208	67	CM-812	69
CM-816	63	GM-34-M	67	CM-924	69
CM-543	63	GM-181	67	CM-954	69
CM-817	63	GMJ33M	67	CM-11-M	69
CM-544	63	GM-324	67	CM-845	69
CM-991	63	GM-183	67	CM-813	69
CM-932	63	GM-281	67	CM-925	69
CM-528	63	GM-291	67	CM-962	70
CM-222	64	GM-282	67	CM-801	70
SM-3625	64	GM-121	68	CM-914	70
SM-3626	64	GM-321	68	CM-964	70
SM-5688	64	GM-69-M	68	CM-966	70
SM-5687	64	GM-122	68	CM-765-L	70
SM-5690	64	GM-67-M	68	CM-767-L	70
SM-5624	64	GM-322	68	CM-903	70
SM-5627	64	GM-123	68	CM-974	70

CM-904	70	8886661231	78	8886612701	83
CM-31-M	70	8886660841	78	8886613911	83
CM-975	70	8886663941	78	8886612711	83
CM-905	70	8886661241	78	8886611511	83
CM-790L	70	8886661831	79	8886614811	83
CM-776L	70	8886661841	79	8886614821	83
CM-779L	70	8886664141	79	8886612111	83
CM-780L	70	8886664151	79	8886614711	83
UM-879	71	8886661851	79	8886612411	83
UM-978	71	8886662651	79	8886616523	83
UM-878	71	8886605941	79	8886614721	83
UM-877	71	8886674441	79	8886612421	83
UM-977	71	CMM-435-MG	79	8886612821	83
CM-558	71	8886674551	79	8886605731	83
UM-345	71	8886674451	79	8886646221	83
UM-245	71	88866751-71	79	8886646231	83
UM-243	71	8886675551	80	8886600111	84
UM-383	71	8886675561	80	8886621311	84
CM-90-M	72	8886675571	80	8886621321	84
CM-436	72	88866757-71	80	8886620133	84
CM-827	73	8886678771	80	8886621331	84
CM-740L	73	8886676971	80	8886613211	84
GMT-226	74	8886653311	81	8886612221	84
CMT-811	74	8886653211	81	8886620721	84
CMT-740-M	74	8886653221	81	8886620731	84
CMT-812	74	SMM-5516	81	8886620741	84
CMT-967	74	SMM-5526	81	8886614921	84
CMT-741-M	74	SMM-5536	81	8886613431	84
CM-951	74	SMM-5221	81	8886621721	84
CM-949-M	74	SMM-5431	81	8886621933	84
CMT-950-M	75	SMM-5042	81	8886621731	84
CMT-974	75	SMM-5043	81	8886621943	84
CMT-949-M	75	SMM-5033	81	8886621741	84
CMT-789-L	75	SMM-5323	81	8886622131	84
		SMM-5044	81	8886603641	84
<b>Maxon™</b>		8886663241	81	8886622933	85
8886600741	78	SMM-5241	81	8886603433	85
8886600751	78	8886662451	81	8886623331	85
8886601563	78	88866704-23	82	8886622943	85
8886660521	78	SMM-5831	82	8886623341	85
8886660531	78	8886690831	82	8886623241	85
8886611831	78	8886690841	82	8886603443	85
8886660821	78	8886639053	82	8886622953	85
8886660831	78	8886614501	83	8886623351	85
8886663931	78	8886613901	83	8886623131	85

8886623531	85	GMM-343L	87	SP-661G	96
8886614331	85	GMM-344L	87	SP-683	96
8886623541	85	8886640351	87	SP-683G	96
8886624541	85	8886640451	87	SP-665G	96
8886624753	85	8886623741	88	SP-665	96
8886624551	85	8886623761	88	SP-655	96
8886624561	85	8886624961	88	SP-629	96
8886624571	85	8886633961	88	SP-684	96
8886626131	85	88866339-71	88	SP-684G	96
8886626141	85	8886645561	88	SP-868G	97
8886626151	85	8886645661	88	SP-868	97
8886626161	85	8886641971	88	SP-685	97
8886626741	86	88866768-11	89	SP-670	97
8886627553	86	88866297-11	89	SP-675	97
8886626751	86	88866297-21	89	SP-624-BC	97
8886626951	86	88866288-11	89	SP-649	97
GMM-951	86	88866288-21	89	SP-689	97
8886627563	86	VMM-137	89	SP-690	97
8886626761	86	88866292-51	89	CP-460	98
8886626961	86	8886627761	89	3-CP-360	98
8886627573	86	GMM-873L	90	CP-443	98
8886626771	86	8886629571	90	CP-453	98
8886626971	86	GMM-806L	90	SP-1696G	99
8886625741	86	GMM-807L	90	SP-1696	99
8886625751	86	GMMT-951	90	SP-1697	99
8886625761	86	GMMT-563-MG	90	SP-1606	99
8886628751	86	GMMT-961	90	SP-5697G	99
8886628561	86	GMMT-541-MG	91	SP-2698G	99
8886628761	86	GMMT-540-MG	91	SP-5695-G	99
8886628571	86	GMMT-89L	91	SP-5695	99
8886628771	86	GMMT-90L	91	SP-5698-G	99
GMM-563L	86	8886606163	91	SP-5698	99
GMM573L	86	8886606173	91	SP-5699G	99
8886631163	87	8886606773	91	SP-5699	99
8886630761	87	GMM-496L	92	SP-5680G	99
8886631173	87	GMM-497L	92	SP-5681G	99
8886630771	87	8886662801	92	SP-5681	99
8886630972	87			SP-5692G	99
GMM-634L	87	<b>Surgipro™ /</b>		SP-1693G	99
GMM-632L	87	<b>Surgipro™ II</b>		SP-5686G	99
8886631962	87	SP-660G	96	SP-5686	99
8886632271	87	SP-660	96	SP-5682G	99
8886631972	87	SP-559G	96	SP-5682	99
GMM-341L	87	SP-631	96	SP-5687	99
		SP-661	96		

SP-5687G	99	VPF-735-X	102	VP-706-X	104
SP-5587	99	VP-738-X	103	VP-718-MX	104
SP-5663G	99	VP-738-MX	103	VP-718-X	104
SP-5563	99	VP-630-X	103	VP-726-MX	104
SP-1617	100	VP-631-X	103	VP-726-X	104
SP-1618	100	VP-701-X	103	VP-76-MX	104
SP-1636	100	VP-701-MX	103	VP-76-X	104
SP-1635	100	VP-702-X	103	VP-78-X	105
SP-1634	100	VP-702-MX	103	VP-705-MX	105
SP-1644	100	VP-713-X	103	VP-705-X	105
SP-1643	100	VP-713-MX	103	VP-70-X	105
SP-621	100	VP-746-X	103	VP-719-X	105
SP-622	100	VP-746-MX	103	VP-720-X	105
SP-622-BC	100	VP-703-X	103	VP-725-MX	105
SP-623-BC	100	VP-703-MX	103	VP-725-X	105
SP-623	100	VP-709-MX	103	VP-75-X	105
SP-8656G	100	VP-709-X	103	VP-719-MX	105
P-2790G	101	VP-365-X	103	VP-897-M	105
P-2794G	101	VP-72-X	103	VP-720-MX	105
P-2795G	101	VP-72-MX	103	VP-813-X	105
VP-900-X	101	VP-37-MX	103	VPF-726-X	105
VP-901-MX	101	VP-73-X	103	VP-126-MX	105
VP-900-MX	101	VP-73-MX	103	VPF-706-X	105
VP-902-X	101	VP-102-MX	103	VPF-720-X	105
VP-903-MX	101	VPF-702-X	103	VP-120-MX	105
VP-902-MX	101	VPF-713-X	103	VP-86-X	105
VP-973-X	101	VPF-703-X	103	VP-881-X	105
VP-904-X	101	VPF-709-X	103	VP-570X	105
VP-904-MX	101	VPF-703-MX	103	VP-580-X	105
VP-209MX	102	VP-763-X	104	VP-584-X	105
VP-744-X	102	VP-875-MX	104	VP-582-X	105
VP-744-MX	102	VP-875-X	104	VP-581-X	105
VP-745-X	102	VP-759-MX	104	VP-583-X	105
VP-743-MX	102	VP-74-X	104	VP-585-X	105
VP-753-MX	102	VP-776-X	104	VP-841-X	105
VP-736-M	102	VPF-776-X	104	VP-851-X	105
VP-628	102	VP-652-MX	104	VP-532-X	105
VP-734-MX	102	VP-889-X	104	VP-533	105
VP-735-X	102	VP-890-X	104	VP-772-X	106
VP-737-MX	102	VP-812-X	104	VP-754-X	106
VP-735-MX	102	VP-727-X	104	VP-755-X	106
VP-747-X	102	VP-800-X	104	VP-854-X	106
VP-747-MX	102	VP-363X	104	VPF-610-X	106
VPF-738-MX	102	VP-706-MX	104	VP-714-X	106

VP-711-X	106	VP-522-X	108	VPF-704-X	111
VP-733-MX	106	VP-422X	108	VPF-807-X	111
VP-733-X	106	VP-542-X	108	VP-911-X	111
VP-715-X	106	VP-523	108	VP-707-X	111
VP-205-X	106	VP-543	108	VP-728-X	111
VP-710-MX	106	VPF354X	108	VP-728-MX	111
VP-710-X	106	VPF-521-X	108	VP-708-X	111
VP-204-X	106	VPF-522-X	108	VP-751-X	111
VPF-733-X	106	VPF-354-X	108	VPF-728-X	111
VPF-710-X	106	VPF-543	108	VP-941-X	112
VPF-711-X	106	CP-16-M	108	VP-942-X	112
VP-972X	106	CP-411	108	VP-914-X	112
VP-870-X	106	CP-17-M	108	VP-924-X	112
VP-871-X	106	CP-412	108	VP-925-X	112
VP-872-X	106	VP-945	109	VP-926-X	112
VP-555X	107	CP-422	109	VP-965	112
VP-660-X	107	CP-423	109	VP-956	112
VP-556X	107	CP-10-M	109	VP-934-X	112
VP-566X	107	CP-424	109	VP-935-X	112
VP-556MX	107	CP-425	109	VP-435-X	112
VPP-557-X	107	VP-842-X	109	VP-936-X	112
VP-557MX	107	VP-843	109	VP-937	112
VP-757-X	107	VP-853	109	VP-975-X	112
VP-557X	107	CP-433	109	VP-976-X	112
VP-558X	107	CP-434	109	VP-986-X	112
VP-559-X	107	CP-435	109	VP-977	112
VPF-557-X	107	CP-535	109	VP-978	112
VP-331X	107	CP-30-M	109	CP-844	112
VP-771-X	107	CP-454	109	CP-845	112
VP-761-X	107	CP-455	109	CP-957	113
VP-761-MX	107	3CP-444	109	CP-836	113
VP-760-X	107	CP-824	109	P-2771K	113
VP-762-X	107	CP-825	109	P-2756K	113
VP-831X	108	CP-824-L	110		
VP-832X	108	CP-834-L	110		
VP-833	108	GP-282	110		
VP-834	108	CP-92-M	110		
VP-510X	108	CP-414	110		
VP-520X	108	CP-415	110		
VP-511-X	108	VP-721-MX	111		
VP-621-X	108	VP-704-X	111		
VP-521-X	108	VP-704-MX	111		
VP-421X	108	VP-712-X	111		
VP-541X	108	VP-712-MX	111		
				<b>Ti-Cron™ / Surgidac™</b>	
				88863008-21	116
				88863008-31	116
				88863003-32	116
				88863008-41	116
				88863001-41	116
				88863003-42	116
				88863011-51	116
				88863008-51	116
				88863003-52	116

88863001-62	116	88863028-81	118	88863226-41	122
88863002-62	116	88863155-81	119	88863186-21	122
88863001-72	116	88863026-71	119	88863186-31	122
88863380-42	116	88863026-81	119	88863186-41	122
88863380-52	116	88863027-79	119	88863186-51	122
88863380-62	116	SCD-3058G	119	88863271-21	122
88863380-72	116	SCD-2864G	119	88863087-31	122
88863380-82	116	SCD-3059G	119	88863271-31	122
88863381-09	116	88862862-34	119	88863084-46	122
88863381-89	116	SCD-3049G	119	88863087-41	122
88863045-21	117	88862854-33	120	88863271-41	122
88863046-31	117	88863109-43	120	88863086-56	122
88863047-31	117	88863059-53	120	88863083-56	122
88863047-41	117	88863127-79	120	88863084-56	122
88863055-41	117	88863163-51	120	88863091-56	122
88863048-51	117	88863163-61	120	88863087-51	122
88863069-41	117	88863163-71	120	88863273-56	122
88863050-41	117	88863062-41	120	88863269-56	122
88863050-51	117	88863062-51	120	88863267-42	123
88863065-51	117	88863165-51	120	88863348-31	123
88863050-61	117	88863062-61	120	88863369-31	123
88863065-61	117	88862947-53	120	88863227-41	123
88863070-51	117	88863288-01	121	88863369-41	123
88863023-51	117	88863106-31	121	88863385-56	123
88863023-61	117	88863280-11	121	88863227-51	123
88863012-61	117	88863314-21	121	88863348-51	123
88863023-71	117	88863280-21	121	88863323-56	123
88863054-51	118	88862801-21	121	88862978-51	123
88863054-61	118	88863314-31	121	88863277-56	123
88863054-71	118	88863280-31	121	88863344-56	123
88863054-81	118	88863280-41	121	88862839-51	123
88863056-89	118	88862636-56	121	88863369-51	123
88863128-79	118	88862775-21	121	88863257-62	124
88863128-82	118	88862775-31	121	88863370-51	124
CD-496	118	88862775-41	121	88863185-31	125
88863395-51	118	88862627-41	121	88863185-41	125
88863395-61	118	88863230-56	121	88863353-56	125
88863395-71	118	88863298-21	122	88863033-52	125
88863395-81	118	88863297-21	122	88863185-61	125
88863309-51	118	88863139-21	122	88863185-71	125
88863309-61	118	88863226-21	122	88863185-51	125
88863309-71	118	88863226-31	122	88863256-31	125
88863154-09	118	88863249-31	122	88863336-41	125
88863028-71	118	88863295-41	122	88863256-41	125

88863035-41	125	88863097-61	128	88863204-51	131
88863324-56	125	88863261-72	128	88863160-51	131
VCD-346	125	88863097-71	128	88862851-56	131
88863356-56	125	88863261-82	128	88863303-56	131
88863040-56	125	88863190-61	128	88863241-56	131
88863037-56	125	88862768-51	128	88863157-51	131
88863336-51	125	88863017-51	128	88863088-41	132
88863256-51	125	88863017-61	128	88863088-51	132
88862981-56	125	88863103-71	128	88863167-31	132
88862914-51	125	88863374-52	129	88863088-41	132
88862953-51	125	88863374-62	129	88863167-41	132
88863307-56	126	88863375-71	129	88863222-56	132
88863306-56	126	88863386-62	129	88863213-56	132
88863220-56	126	88863379-52	129	88863203-56	132
88863035-51	126	88863372-72	129	88863202-51	132
88863256-61	126	88863372-82	129	88863088-51	132
88862914-61	126	88863113-71	129	88862855-51	132
88863258-32	126	88863113-81	129	88863290-56	132
88863258-42	126	88863111-79	129	88863212-51	132
88863258-52	126	88863376-81	130	88863167-51	132
88863129-51	126	88863107-31	130	88862815-51	132
88863015-61	126	88863182-31	130	88862852-56	132
88863251-51	126	88863182-41	130	88862840-51	132
88863093-51	126	88863365-56	130	88863213-66	132
88863259-62	126	88863233-51	130	88863202-61	132
88863093-61	126	88863211-56	130	88863167-61	132
88863205-41	127	88863182-51	130	88863131-56	133
88863205-51	127	88863052-56	130	88863275-56	133
88863205-61	127	88862658-56	130	88863159-31	133
88863082-61	127	88862631-51	130	88863159-41	133
88863082-51	127	88863218-36	131	88863236-56	133
88863274-51	127	88863160-31	131	88863235-56	133
88863393-51	127	88862880-41	131	88862879-56	133
88863274-61	127	88863160-41	131	88863159-51	133
88863393-61	127	88863221-56	131	88863193-51	133
88863260-52	127	88863219-56	131	88863319-51	133
88863090-51	127	88863218-56	131	88863240-56	133
88863260-62	127	88863125-51	131	88863228-56	133
88863090-61	127	88862880-51	131	88863239-51	133
88863260-72	127	88862877-56	131	88863168-51	133
88863090-71	127	88862811-51	131	88863119-61	133
88863146-81	127	88862808-51	131	88863092-61	134
88863090-81	127	88863191-51	131	88863092-71	134
88863261-62	128	88863255-56	131	88863147-83	134



88863147-89	134	88864750-31V	139	88862382-83	148
88862818-89	134	88864730-51V	139	88862466-25	148
CD-3113K	135	88864731-31V	139	88862414-89	148
CD-3123K	135	88864731-41V	139	88862393-89	148
D-1764K	135	88864731-51V	139	88862376-89	148
D-1763K	135	88864734-11V	139	88862392-69	148
D-1722K	135	88864734-21V	139	88862393-69	148
D-1760K	135	88864714-11V	139	88862392-49	148
D-1683K	135			88862396-89	149
CD-0423K	135	<b>Flexon™</b>		88862497-89	149
D-1780K	135	88862597-63	142	88862495-89	149
D-1779K	135	88862623-53	142	88862494-89	149
		88862617-53	142		
<b>Vascufil™</b>		88862615-53	142	<b>Dermalon™ /</b>	
88864741-21V	138	88862589-63	142	<b>Monosof™</b>	
88864788-05V	138	88862592-43	143	N-63	152
88864795-06V	138	88862591-43	143	N-59	152
88864703-01V	138	88862586-43	143	N-64	152
88864705-01V	138	88862593-63	143	SN-770G	152
88864702-05V	138			SN-667G	152
88864704-01V	138	<b>Steel</b>		SN-667	152
88864701-05V	138	88862407-53	146	88861741-11	152
88864703-11V	138	88862407-63	146	SN-367	152
88864705-11V	138	88862408-73	146	SN-771G	152
88864704-11V	138	88862408-83	146	SN-871	152
88864756-11V	138	88862407-43	146	88861741-21	152
88864722-21V	138	88862408-73	146	SN-368	152
88864707-11V	138	88862413-83	146	SN-772G	152
88864709-11V	138	88862420-83	146	SN-872	152
88864710-11V	138	88862420-09	146	88861741-31	152
88864709-21V	138	88862410-89	146	88861747-11	153
88864710-21V	138	88862424-89	146	SN-660	153
88864711-21V	138	88862412-89	146	SN-760	153
88864739-21V	138	88862410-69	146	SN-659	153
88864739-31V	138	88862224-89	147	SN-659G	153
88864716-21V	138	88862222-89	147	SN-761	153
88864716-31V	138	88862226-89	147	SN-762	153
88864742-15V	139	88862224-69	147	SN-763	153
88864741-11V	139	88862226-69	147	88861756-11	153
88864794-01V	139	88862222-69	147	88861756-21	153
88864719-21V	139	88862224-49	147	SN-661G	153
88864719-31V	139	88862228-89	147	SN-661	153
88864750-51V	139	88862345-53	148	SN-643	153
88864750-41V	139	88862348-53	148	SN-552G	153
				SN-662G	153

SN-662	153	SN-673G	155	SN-5668	157
88861756-31	153	SN-673	155	SN-5694G	157
88861757-31	153	88861727-41	155	SN-5694	157
SN-644	153	88861727-51	155	SN-5665G	158
SN-653	153	SN-674	155	SN-5666G	158
88861756-41	153	SN-775	155	SN-5666	158
SN-641	153	SN-776	155	SN-5677	158
88861757-41	153	SN-777	155	SN-5667G	158
SN-665	153	CN-290	155	SN-5667	158
88861756-51	153	CN-490	155	SDN-5666G	158
SN-750	153	SN-629	156	SN-5669G	158
SN-666	153	SN-630	156	SN-5669	158
SN-658	154	SN-631	156	SN-5670	158
88861744-31	154	SN-247	156	SN-5862G	159
88861799-31	154	SN-681	156	SN-5663G	159
SN-649	154	SN-781	156	SN-5663	159
SN-663G	154	SN-682	156	SDN-5943G	159
SN-663	154	SN-782	156	SN-5660G	159
88861744-41	154	SN-881	156	SN-5662G	159
88861799-41	154	CN-819	156	SN-5662	159
SN-669	154	CN-822	156	SN-1956	160
SN-648	154	SN-1647G	157	SN-1955	160
SN-758	154	SN-1647	157	SN-1954	160
88861799-51	154	SN-3686	157	SN-1966	160
SN-747	154	SN-3697	157	SN-3965	160
SN-764G	154	SN-1689G	157	SN-1964	160
SN-764	154	SN-5697	157	SN-1995	160
88861744-51	154	SN-1689	157	SN-1994	160
SN-664G	154	88861754-13	157	SN-1693	160
SN-664	154	SN-3695	157	88861742-21	160
88861720-51	154	SN-5696	157	88861742-31	160
SN-765	154	SN-5696G	157	SN-568	160
SN-670	154	SDN-5696G	157	88861742-41	160
88861750-41	154	SN-5698G	157	CN-791	160
SN-671G	154	SN-5690G	157	CN-792	160
SN-671	154	SN-5698	157	CN-793	160
SN-771	154	SN-5690	157	SN-615	161
SN-688	154	88861749-24	157	SN-616	161
SN-788	154	SN-5699G	157	SN-617	161
SN-689	154	SN-5699	157	SN-627	161
SN-790	154	SN-5691	157	88861764-41	161
SN-693	155	SDN-5691G	157	SN-628	161
SN-694	155	SN-5767G	157	88861764-51	161
SN-672	155	SN-5668G	157	SN-728	161

N-2500	162	N-2760K	167	8886440013	172
N-2510	162	N-2766K	167	8886440023	172
N-2512	162	N-2786K	167	SPB-5413G	172
N-2540	162	N-2787K	167	SPB-5142G	172
N-1603G	162	N-2719K	167	8886402821	172
N-2541	162	N-2716K	167	SPB-5143G	172
N-2514	162	N-2966K	167	8886441431	172
N-2530	162	N-2967	167	8886402831	172
N-2532	162	N-2761K	167	8886440113	172
N-2544	162			8886440123	172
N-2546	162	<b>Novafil™</b>		8886440133	172
N-2543	162	8886441913	170	SPB-1623G	173
N-2548	162	8886441923	170	SPB-5633G	173
N-2547	162	8886441933	170	8886442531	173
CN-225	163	88864555-23	170	SPB-1643G	173
CN-723	163	8886442013	170	8886441543	173
CN-722	163	8886442023	170	8886440223	173
CN-825-L	163	8886442033	170	8886440233	173
CN-724	164	8886403231	170	8886440243	173
CN-830-L	164	8886442043	170	SPB-5833G	173
CN-824	164	88864554-41	170	SPB-5843G	173
CN-825	164	8886442233	170	88864582-43	173
CN-823-L	164	8886442243	170	88864418-53	173
GN-283	164	8886400241	170	88864582-53	173
GN-284	164	8886442253	170	8886440333	173
CN-926	165	8886442153	170	8886440343	173
N-2513	165	8886442431	171	SPB-5433G	173
N-2531	165	8886442441	171	SPB-1945G	174
N-2533	165	8886442451	171	88864490-41	174
CN-771	165	8886442461	171	88864490-51	174
N-2714K	166	8886442641	171	XNF-697	174
N-2781K	166	8886442651	171	XNF-1966	174
N-2730-K	166	8886442661	171	XNF-1965	174
N-2770-K	166	8886400963	171	XNF-1964	174
N-2771-K	166	8886443371	171	XNF-1995	174
N-2736K	166	8886443671	171	XNF-1994	174
N-2718K	166	CPB-738G	171	XNF-1993	174
N-2717-K	166	8886441003	172	8886445241	175
N-2799-K	166	8886441013	172	8886445251	175
N-2705K	166	8886441023	172	8886445261	175
N-2704K	166	SPB-1213G	172	8886445951	175
N-2734K	167	SPB-5223G	172	8886445853	175
N-2757K	167	8886402721	172	8886445863	175
N-2756K	167	SPB-1233G	172	8886445961	175

8886445971	175	88861935-61	181	88861964-61	184
88864619-51	175	88861991-42	181	88861962-52	185
8886445041	175	88861883-51	181	88861971-61	185
8886445463	175	88861883-61	181	88861962-72	185
8886445561	175	88861883-71	181	88861971-71	185
8886445473	175	88861922-71	181	88861998-62	185
8886445571	175	88861932-51	181	88862011-51	185
88864463-71	175	88861896-71	182	88861928-52	185
88864465-81	176	88861915-61	182	88861928-62	185
88864464-73	176	88861915-71	182	88861928-72	185
88864465-71	176	88861913-72	182	88861909-62	185
88864538-63	176	88861895-71	182		
88864470-81	176	SBS-1928G	182	<b>Sofsilik™</b>	
8886440861	176	SBS-1880G	182	S-182	188
8886446063	177	88861981-33	182	S-193	188
PB-7740K	177	SBS-1884G	182	S-183	188
PB-6713K	177	SBS-1883G	182	S-243	188
PB-6723K	177	88861985-51	183	S-303	188
88861918-31	180	88861985-61	183	S-194	188
88861919-31	180	88861985-71	183	S-184	188
88861899-41	180	88861948-61	183	S-254	188
88861918-41	180	88861948-81	183	S-244	188
88861919-41	180	88861910-22	183	S-404	188
88861899-51	180	88861910-32	183	S-304	188
88861918-51	180	88861957-32	183	S-195	188
88861919-51	180	88861927-32	183	S-185	188
88861917-61	180	88861995-31	183	S-245	188
88861902-61	180	88861927-42	183	S-305	188
88861919-61	180	88861995-41	183	S-605	188
88861917-71	180	88861989-32	184	S-196	188
88861902-71	180	88861989-42	184	S-176	188
88861919-71	180	88861993-62	184	S-346	188
88861917-82	180	88861958-42	184	S-246	188
88861906-42	180	88861951-41	184	S-316	188
88861906-52	180	88861926-42	184	S-606	188
88861920-62	180	88861958-52	184	S-197	188
88861920-72	180	88861967-52	184	S-187	188
88861914-83	180	88861951-51	184	S-317	188
88861920-93	180	88861951-61	184	S-607	188
88861931-43	181	88861969-52	184	S-318	188
88861982-51	181	88861969-62	184	S-608	188
88861935-41	181	88861968-61	184	S-282	189
88861934-53	181	88861969-72	184	S-610	189
88861935-51	181	88861964-51	184	S-206	189

S-199	189	SS-675	192	SS-722	196
LS-640	189	SS-677	192	SS-623	196
LS-639	189	SS-679	192	SS-723	196
LS-638	189	SS-678	192	SS-624	196
LS-637	189	SS-70-M	192	VS-802	196
LS-636	189	SS-680	192	VS-809	196
S-1768K	189	SS-689	192	VS-889	197
S-1769K	189	SS-648	193	VS-890	197
S-1746K	189	SS-647	193	VS-891	197
S-1750K	190	SS-646	193	VS-806	197
S-1172	190	SS-651	193	VS-706-M	197
S-1173	190	CS-482	193	VS-709	197
S-1174	190	CS-485	193	VS-880	197
S-1373	190	S-1766K	193	VS-881	197
S-1272	190	S-1765K	193	VS-882	197
S-1735	190	SS-82-M	193	VS-823	197
S-1274	190	SS-733	194	VS-581	197
S-1789K	190	SS-732	194	CS-210	198
S-1783K	190	SS-632	194	CS-211	198
SS-681	191	SS-734	194	GS-44-M	198
SS-673	191	SS-633	194	VS-870	199
SS-686	191	CS-390	194	VS-871	199
SS-645	191	CS-490	194	GS-47-M	199
SS-745	191	SS-5639	195	VS-872	199
SS-746	191	SS-1639G	195	VS-873	199
SS-682	191	SS-5649G	195	GS-451	199
SS-652	191	SS-5640	195	GS-452	199
SS-683G	191	SS-5641G	195	GS-453	199
SS-683	191	SS-5641	195	GS-34-M	199
SS-653	191	SS-5678	195	GS-33-M	199
SS-694	191	SS-5676	195	GS-30-M	199
SS-654	191	SS-5677	195	VS-552	199
SS-655	191	SS-5679	195	CS-562	199
SS-656	191	SS-5684	195	GS-64-M	200
SS-629	192	SS-5685G	195	GS-831	200
SS-684	192	SS-1694G	195	GS-63-M	200
SS-684G	192	CS-791	196	GS-67-M	200
SS-784	192	CS-792	196	GS-832	200
SS-695	192	CS-793	196	GS-62-M	200
SS-785	192	CS-794	196	GS-833	200
SS-786	192	SS-523	196	GS-66-M	200
SS-787	192	SS-525	196	GS-834	200
SS-685G	192	SS-621	196	GS-835	200
SS-685	192	SS-622	196	VS-533	200

CS-16-M	200	VLOCM0005	206	VLOCM0324	209
CS-17-M	200	VLOCM0015	206	VLOCM0305	209
GS-822	201	VLOCM0025	206	VLOCM0345	209
GS-823	201	VLOCM0035	206	VLOCM0315	209
GS-824	201	VLOCM0113	207	VLOCM3225	209
CS-423	201	VLOCM0123	207	VLOCM1704	209
CS-10-M	201	VLOCM0133	207	VLOCM1744	209
CS-424	201	VLOCM0114	207		
CS-425	201	VLOCM0124	207	<b>V-Loc™ 180</b>	
VS-842	201	VLOCM0134	207	VLOCL2815	212
VS-843	201	VLOCM0115	207	VLOCL2826	212
VS-844	201	VLOCM0125	207	VLOCL1413	212
VS-845	201	VLOCM0135	207	VLOCL0003	212
CS-745	201	VLOCM1824	207	VLOCL0013	212
CS-744	201	VLOCM1624	207	VLOCL0023	212
VS-846	201	VLOCM0223	207	VLOCL0033	212
CS-746	201	VLOCM0224	207	VLOCL0004	212
CS-93-M	202	VLOCM1203	208	VLOCL0014	212
CS-92-M	202	VLOCM1204	208	VLOCL0024	212
CS-85-M	202	VLOCM0843	208	VLOCL0034	212
X-1095	202	VLOCM0813	208	VLOCL0015	212
VS-863	202	VLOCM1904	208	VLOCL0025	212
VS-964	202	VLOCM0804	208	VLOCL0123	213
VS-766G	202	VLOCM1944	208	VLOCL0133	213
S-2176K	203	VLOCM0844	208	VLOCL0114	213
S-1780K	203	VLOCM0814	208	VLOCL0124	213
S-2780K	203	VLOCM0824	208	VLOCL0134	213
S-2782K	203	VLOCM0603	208	VLOCL0115	213
S-2752K	203	VLOCM0623	208	VLOCL0125	213
S-1734K	203	VLOCM0604	208	VLOCL0135	213
S-1733K	203	VLOCM2004	208	VLOCL0223	213
S-1732K	203	VLOCM0644	208	VLOCL0224	213
		VLOCM2044	208	VLOCL0803	213
<b>V-Loc™ 90</b>		VLOCM0614	208	VLOCL0813	213
VLOCM1413	206	VLOCM0624	208	VLOCL0804	213
VLOCM1423	206	VLOCM0625	208	VLOCL1904	213
VLOCM0003	206	VLOCM2105	209	VLOCL0844	213
VLOCM0013	206	VLOCM2205	209	VLOCL0814	213
VLOCM0023	206	VLOCM2145	209	VLOCL0824	213
VLOCM0033	206	VLOCM2245	209	VLOCL0603	214
VLOCM0004	206	VLOCM2115	209	VLOCL0613	214
VLOCM0014	206	VLOCM1544	209	VLOCL0604	214
VLOCM0024	206	VLOCM1545	209	VLOCL0644	214
VLOCM0034	206	VLOCM0344	209	VLOCL0614	214

VLOCL0624	214	VCD20305MFE	222
VLOCL0615	214	VCD20305MSE	222
VLOCL0625	214	VCD20305FE	222
VLOCL2105	214	VCD20305SE	222
VLOCL2145	214	VCD2031MFE	222
VLOCL1544	214	VCD2031MSE	222
VLOCL1545	214	VPG36FE	222
VLOCL0344	214	VPG36SE	222
VLOCL0324	214	VPG56FE	222
VLOCL0305	214	VPG56SE	222
VLOCL0345	214	VCD20331FE	222
VLOCL0315	214	VCD2036305FE	222
VLOCL0325	214	VCD2030556MFE	222
VLOCL0335	214	VCD2030556MSE	222
VLOCL0306	214	8886332356	222
VLOCL0346	214	8886332456	222
VLOCL0316	214	8886321956	222
VLOCL0326	214	8886330656	222
VLOCL0336	214	8886309156	222
VLOCL0416	214	8886321356	222
VLOCL0426	214	VCD656	222
VLOCL0436	214	8886308656	222
VLOCL2744	215	8886321856	222
		8886303756	222
<b>V-Loc™ PBT</b>		VCD110	222
VLOCN0814	218	8886867701	222
VLOCN0604	218		
VLOCN0644	218	<b>Wound Closure accessories</b>	
VLOCN0614	218	88868619-03	223
VLOCN0605	218	88868618-00	223
VLOCN0615	218	BW25G	223
VLOCN1126	218	88868631-12	223
VLOCN0305	219		
VLOCN0325	219	<b>Needle chart</b>	
VLOCN0306	219	Cardiovascular	224
VLOCN0326	219	Plastic/Cosmetic	225
VLOCN0327	219	Cuticular	225
VLOCN0437	219	General Surgery	226
		Microsurgery	226
<b>Pledgets</b>		Urology	226
VCD20316MFE	222	Orthopedic	227
VCD20316MSE	222	Multi-purpose	227
VCD20316FE	222	Ophthalmic	227
VCD20316SE	222		

# Medtronic

IMPORTANT: Please refer to the package insert for complete instructions, contraindications, warnings and precautions.

© 2019 Medtronic. All rights reserved. Medtronic, Medtronic logo and Further, Together are trademarks of Medtronic.

<sup>TM\*</sup> Third party brands are trademarks of their respective owners. All other brands are trademarks of a Medtronic company.  
16-eu-or-ess-sutures-catalogue-dic16-897730-fr

[medtronic.com/covidien/en-gb/index.html](https://www.medtronic.com/covidien/en-gb/index.html)