

- 48 Years, China Longest Time APM Factory
- 30,000 sq.m., China Biggest APM Factory Zone
- Japan TOSHIBA, Global Most Professional HF Device



Please check the details with our company profile video:









Taiwan AnnYih Synchronous Motor Alternating Control System



2.8" Bubble Pad Stage I Overlay System





QDC-303 is designed for the treatment and prevention of pressure sore stage I in short-term and domiciliary therapies. The mattress in PVC is composed of 130 7cm-high bubble cells, particularly comfortable, it is fixed to the bed by extra flaps on both top and bottom side. Very easy to install, maintain and use thanks to the possibility of adjust comfort range. The pump can be hung to the end of the bed by means of two hooks.



Without Flaps Optional











20' Container:

Color Box Code: **B2016Y1**


2760 sets (Bubble Pad & Pump)

40' HQ:

6300 sets (Bubble Pad & Pump)




Worldwide Certification

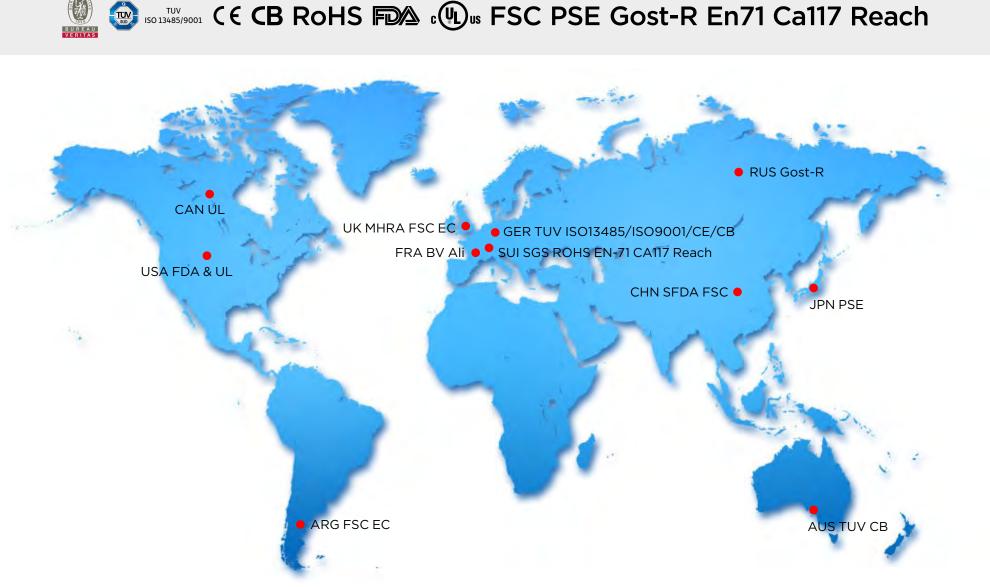





Beige Blue Purple QDC-303 Introduce

0.32mm Taiwan Nanya K80 PVC

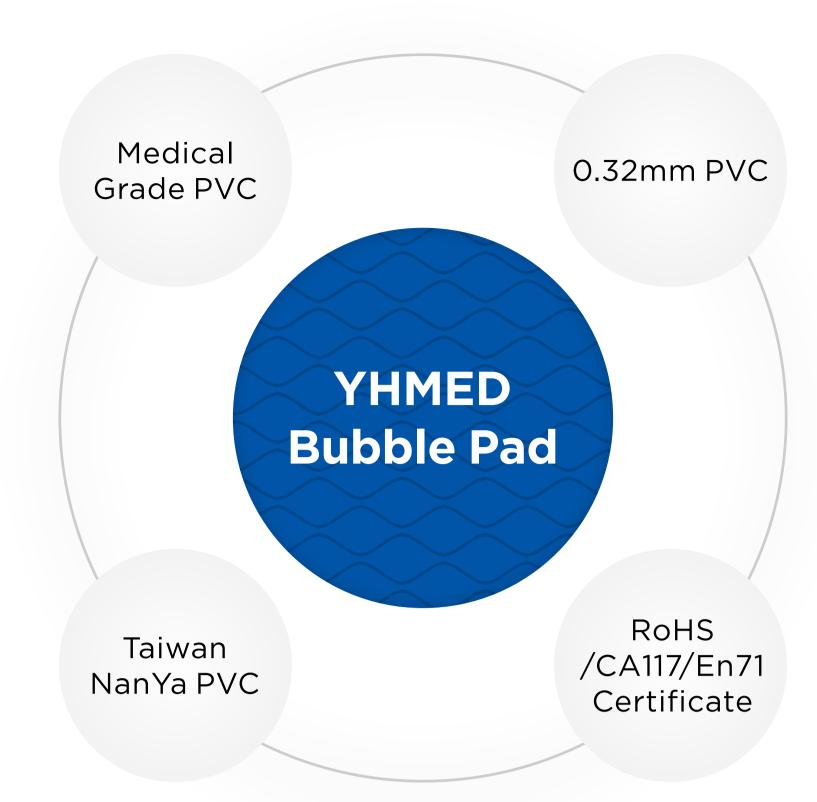

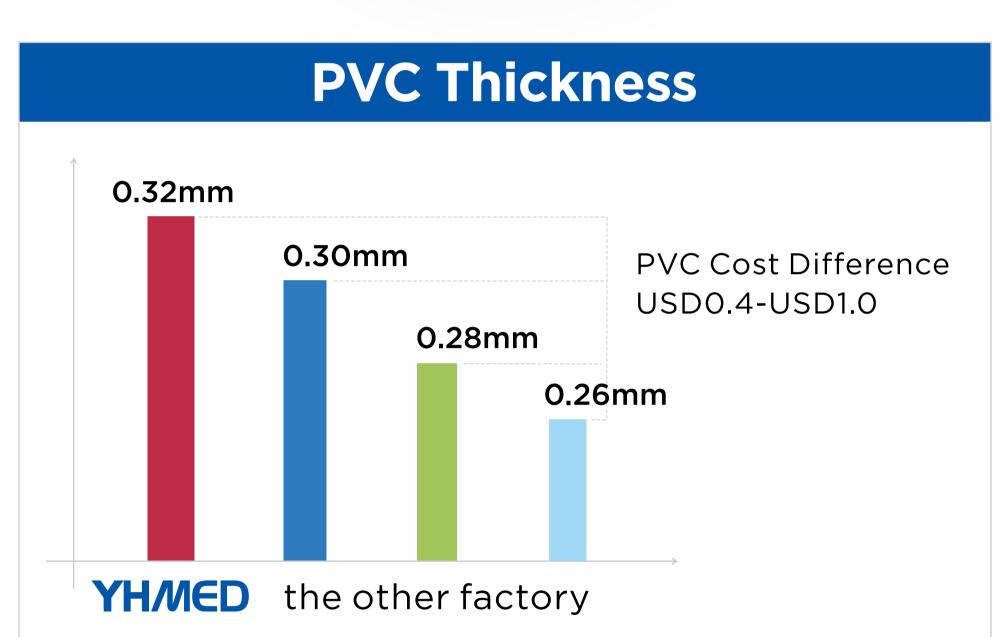

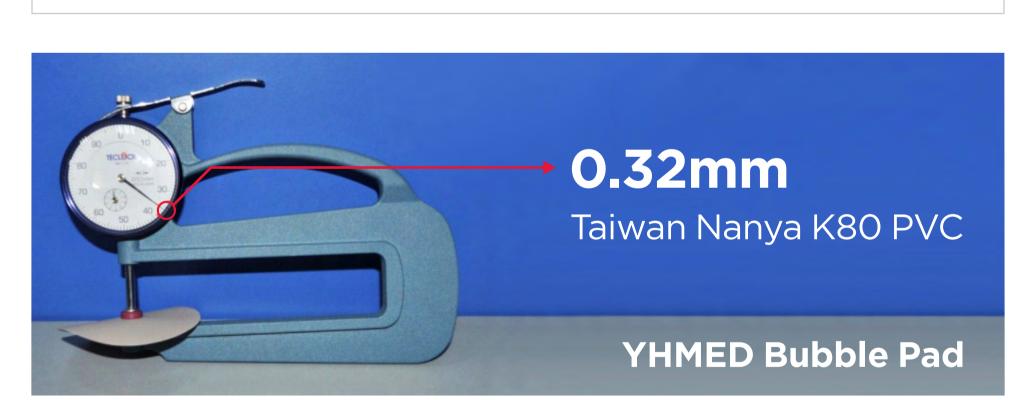

Comparison	YHMED	APEX TW	TW/China PVC Factory
PVC Thickness	0.32 mm	0.28mm-0.3mm	0.23mm-0.28mm
PVC Quality Grade	TW NanYa K80 En-71 Odorless Medical-Grade PVC	TW Nanya Pungent Odor PVC	Pungent Odor PVC





Japan TOSHIBA Power Tube

Global Most Advanced & Biggest 4 Automatic Manufacturing Platforms

200kw One-Time Forming HF Device





Bubble Pad Production Capacity:

1,000,000

Pieces Per Year.



Total Quality Control (TQC)



100% Inflatable Testing Rate,

24 HoursQuality Testing Time.

200×90cm Full Size
Light Box Quality Inspection



50°C High TemperatureLong Term Test



150kg

Rolling Support Weight Test







Free Repairing Kit



YHMED SEER

Technical Parameters Of The Bubble Mattress



Dimensions without flaps	200×90×7 cm (L×W×H)
Extension Flaps length	50 cm+ 50cm (Head & Foot Side)
Mattress material	TW NanYa EN71 K80 Medical-Grade PVC
Material thickness	0.32 mm
Cold-resistamt	-30°C
Number of cells	130
Weight support	135 kg
HF Welding	One-Time Forming (Japan TOSHIBA Power Tube)
Certificate	Rohs / CA-117 / EN-71
Warranty	12 Months



Japan Panasonic Piston Hat












New Economy Compressor
Exciting Competitive Price!



the Most Quiet Pump in China



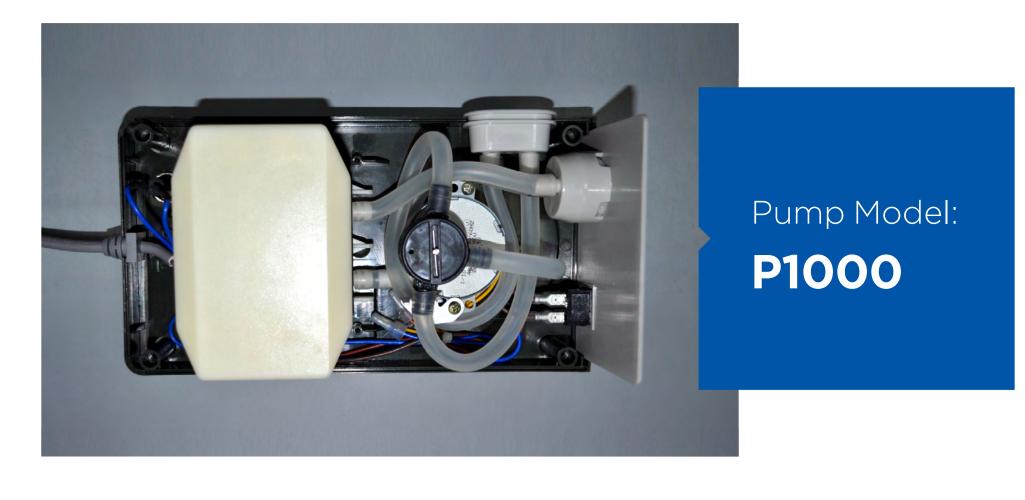

Taiwan AnnYih Synchronous Motor Alternating Control System The Alternating Synchronous Motor that YHMED used is from

Taiwan AnnYih Brand (the repairing rate is less 0.01% within one year).

And the other suppliers that used is the Chinese Synchronous Motor (the repairing rate is 30%-40% within one year).


The cost difference is USD0.5-USD0.6.



Comparison	YH∕M€D	APEX TW	TW/China PVC Factory
Synchronous Motor (Alternating Function)	ANN YIH Global Top-Class Quality Brand	TW Motor	High Failure Rate Stop Alternating
Repair Rate (1 year Used)	0.01%	1~5%	20~50%

The new adding Air Storage for

Japan Panasonic Piston Hat Compressor

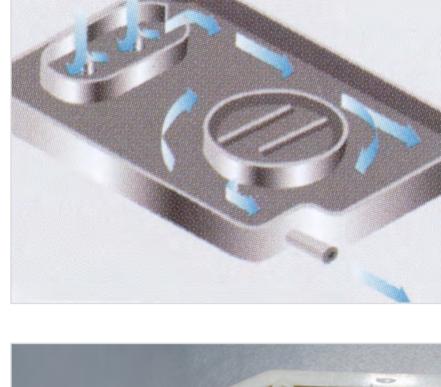

can reduce the working noise of the compressor, more steady output pressure, provide more comfortable for the patient.

YHMED Compressor

Weight:


500g+

Noise-Reduce Technology which




- Japan Panasonic Massage Armchair Piston Hat The cost difference is USD2.5-USD3.0.
- F-400
- 000gX1g/177ozX0 Better than most of China Aluminium Compressor (only 300-400g) The Air Compressor that YHMED used (the compressor used in model P1000) SF-400 is with 2 years Warranty. And the compressor that other suppliers used is with 1 year warranty. RECOMMEND 100g+

CAPACITY:

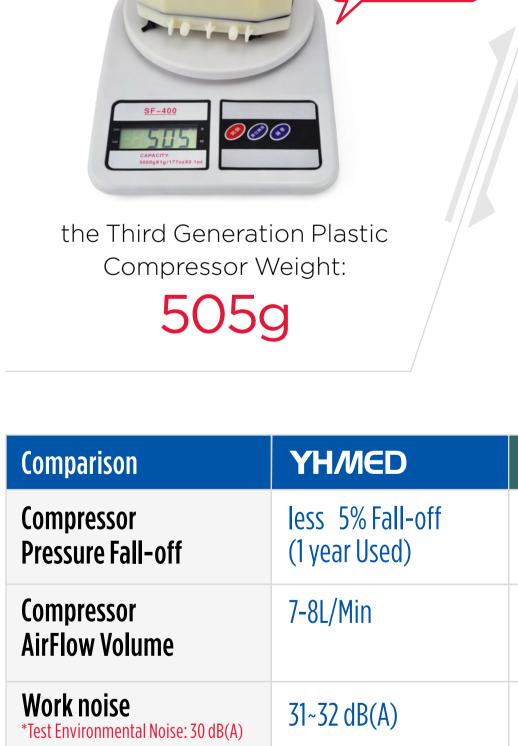

Life Time

(Based on actual test data)

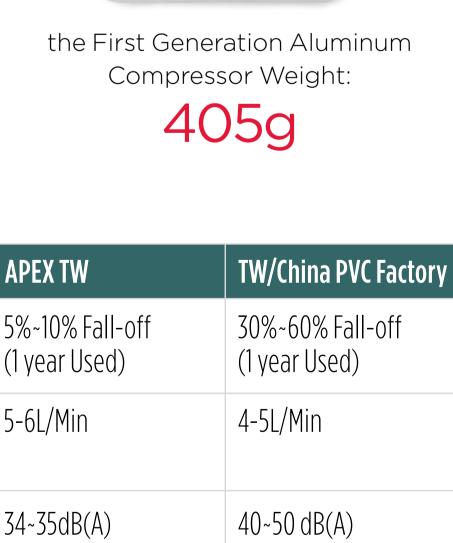

3~12 Months

∅∅


Precise Stable-Comfort-Pressure -Control Technology Better Quality, More Precise & More Stable than PCBA Pressure Control

24 ~ 36 Months



36~48 Months



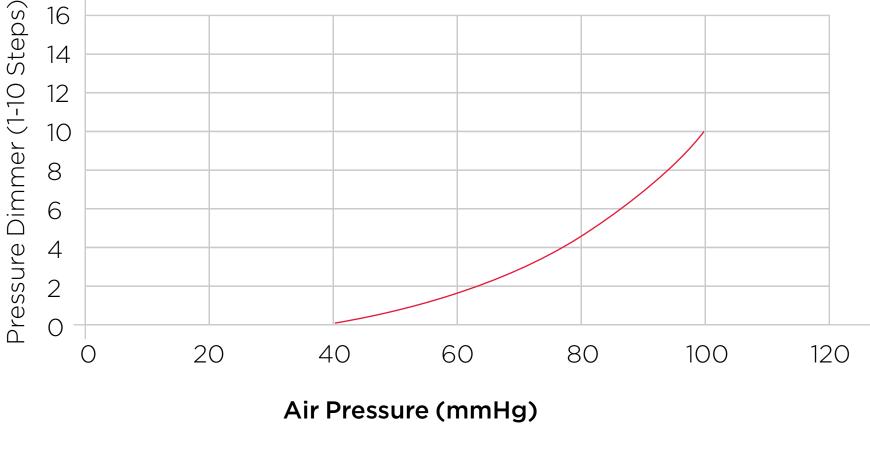



The Actual Data for the Compressor Working Noise

Some Chinese factories always advertise in their catalogue that their compressor (with the first or second generation technology) can reach 18/16dB (A) even just 14dB (A). But this is IMPOSSIBLE! Even for the latest Third Generation Compressor, its noise just can reach 31-32dB (A) which is tested in the authoritative noisy laboratory. If the average noise could reach 31-32dB (A), it will be the excellent lownoise level already.

18 dB(A) is a Joke!



40-50dB(A)Real Sound Level 48dB(A)

- Test Environmental Noise 30 dB(A)

- $= 18 \, dB(A)$
- 10 GB(A)

First Generation Compressor Real Sound Level




YHMED Third Generation Compressor Real Sound Level



* Test Environmental Noise:

30 dB(A)

Almost the Most Quiet PUMP

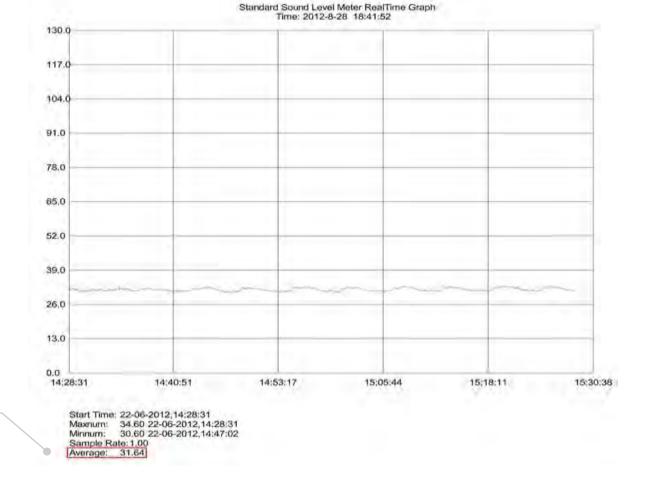



Pump Work noise / Sound Level (Based on actual test data) YH/MED (P1000) China PVC Factory Always show below 20dB(A) Actually 40~50 dB(A) *Test Environmental Noise: 30 dB(A) *Real Sound Level 48dB(A) - Test Environmental Noise 30 dB(A) = 18 dB(A)



Technical Parameters Of The Pump

	P1000	P2000IIE	P4000IIE	
Electrical	AC 110V/220V 50Hz/60Hz			
Pressure range	40-100 mmHg			
Air output	7-8 L/min	6-7 L/min 6-7 L/n		
Dimensions (L×W×H)	26.5×13.5×10 cm	27×13×10.5 cm	25×12×10 cm	
Plastic Casing	Taiwan Chimei ABS			
Test Duration	24 Hours			
Certificate	German TUV-CE IEC/EN60601-1&60601-1-2			
Synchronous Motor	Taiwan Ann Yih Brand (Alternating Control System)			
Alternating	1-in-2			
Warranty	24 Months	12 Months	12 Months	







BLUE IS HOPE

Color Box Package Optional

OEM MOQ: 1000 Sets



29.5×16.5×17.5cm



31×28.5×11cm





AHWED 意志



GREEN IS GROWTH

YHMED E






134819001 (€ CB RoHS PA . . .

EverHealth

Alternating Bubble Mattress

with Adjustable Pump System



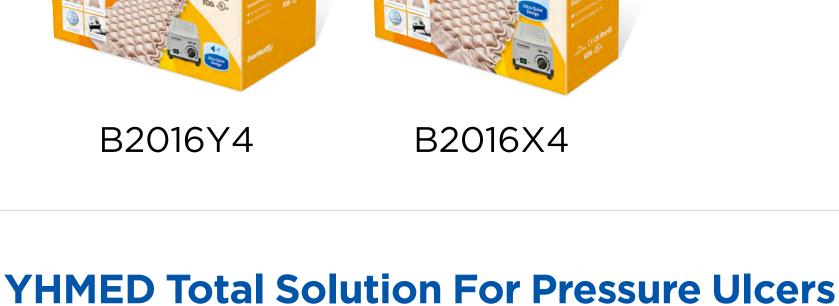

EverHealth.

Alternating Bubble Mattress With Adjustable Pump System

2.8" Bubble

(PVC)



Stage I Stage II Stage III Stage IV

3" /4" /4.5" Cell

(Nylon PVC)

QDC-300B Mattress

+ P4000IIE Pump

3"

Moderate Care Basic Care 5" Nylon PVC 4"

For higher risks of pressure ulcer with

increasing product performance



QDC-303 Beige Mattress

+ P1000 Pump

QDC-303 Beige Mattress + P2000IIE Pump

QDC-303 Blue Mattress

+ P2000IIE Pump


QDC-303 Beige Mattress

+ P4000IIE Pump

QDC-303 Blue Mattress

+ P4000IIE Pump

Features / Benefits

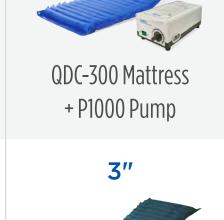

QDC-500 Mattress + P1000 Pump 3"



4.5"

QDC-301 Mattress

+ P1000 Pump

+ P1000 Pump



5" TPU

QDC-5050 Mattress

+ P3000A2QB Pump

5" Nylon TPU

5" / 5" + 3" Cell

QDC-5010E Mattress

+ P3000N2EB Pump

5" Nylon PVC

QDC-5010 Mattress

+ P3000A2QB Pump

Intensive Care





QDC-5080 Mattress

+ P5000 Pump

QDC-5050 Mattress

+ P5000 Pump



QDC-8010 Mattress

+ P3000A2QB3 Pump

TPU

8" Cell



+ P3000A2QB3 Pump

Nylon TPU





Nylon TPU

QDC-8010 Mattress

+ P5000 Pump

TPU



•	
•	
•	
•	
•	

1 40001121	чир	1 10001 dilip	1 3000 1 01	Ч	1 3000 Fallip
Water Matt Seat Cush	_	4.5" C-501B Mattress P2000III Pump			
Air Pump Optional Fo		S Overv	iew		
Air Pump	P1000 Series Size: 26.5×13.5×10cm	P2000 Series Size: 27×13×10.5cm	P3000 Series Size: 30×13.5×11cm	P4000 Series Size: 25×12×10cm	P5000 Series Size: 32×13×21cm
1-in-2 Alternate Mode	•	•	•	•	•
1-in-3 Alternate Mode		•	•		
Static Mode	•		•	•	•
Audio Low Pressure Alarm			•		•
Visible Low Pressure Alarm			•		•
Max Firm Mode					•
YH Quick Connector			•		•
CPC Quick Connector			•		•
Seat Inflate Mode					•
Power Failure Alarm					•
Tech Support Alarm					•
Alarm Mute					•

YHMED 粤華

Since 1969

 48 Years, China Longest Time APM Factory • 30,000 sq.m., China Biggest APM Factory Zone





YEARS

Since 1969

President: Su ZhaoSheng Founded in: 1969

Factory: 30,000 m² 200

Headquarter: Shantou, Guangdong, China

Employees: Sales Network: Over 60 Countries **Production facility:**


Over 100 sets 10-200kw Japan TOSHIBA Power Tube HF Device Manufacturing capacity: Over 1 million sets alternating pressure mattress per year

Heading


Year

APM Milestones We are the leading manufacturer in China with experience in making air mattresses more than 48 years

History & Accomplishments

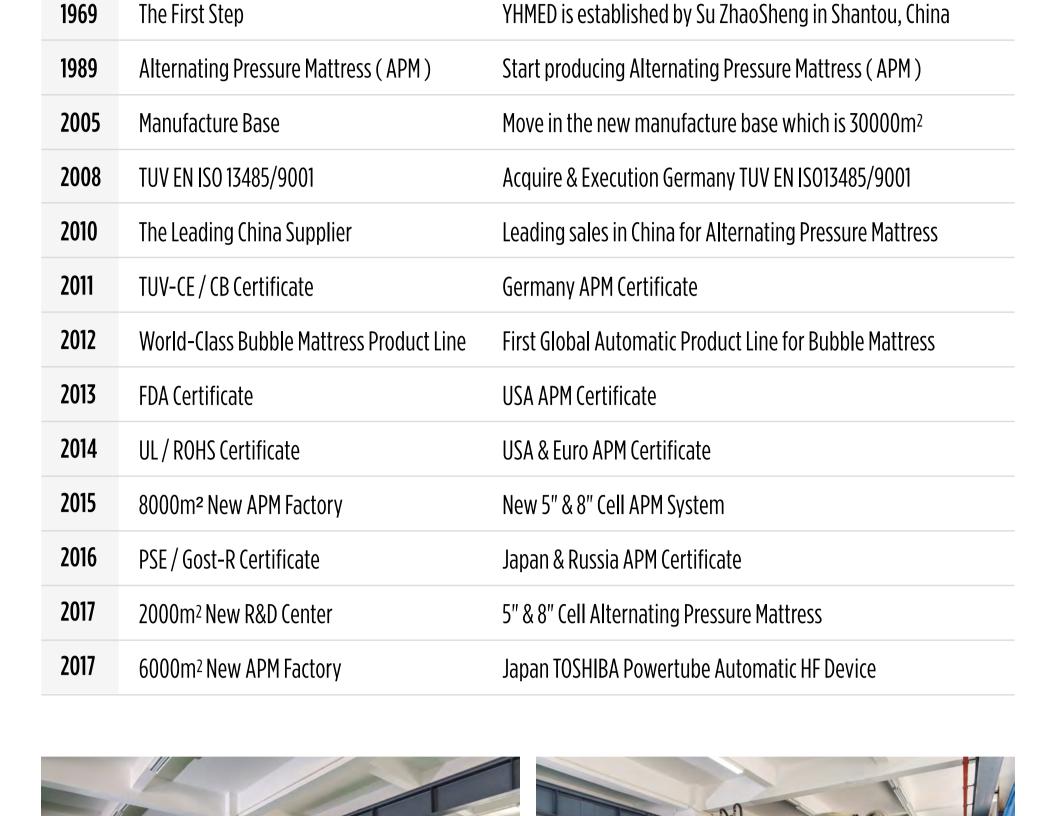




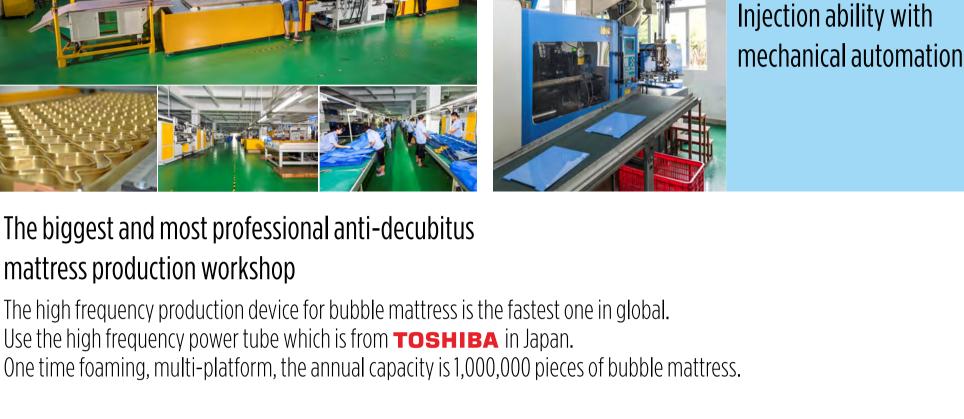

Injection ability with high precision digital control



The pump assemble workshop which is with the most professional automatic line in China.

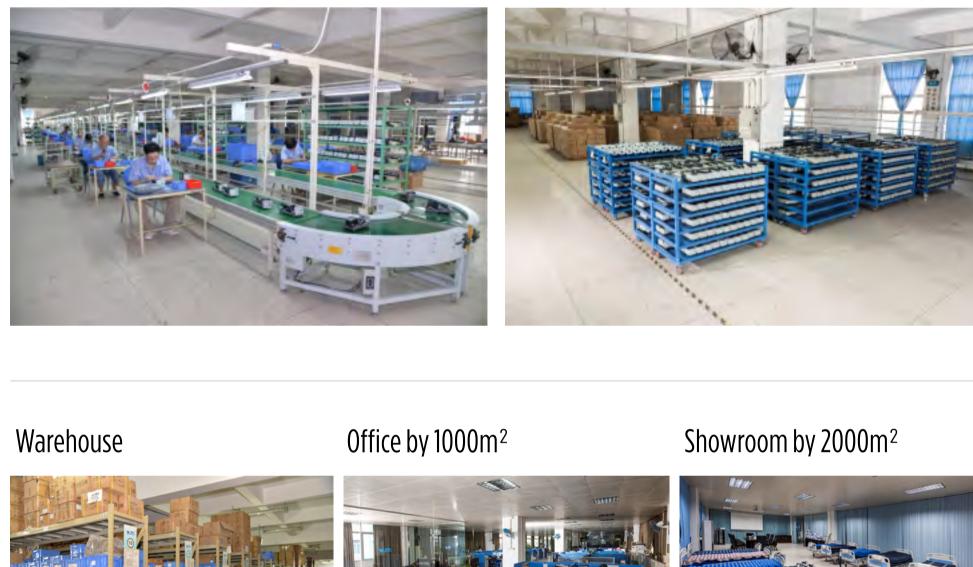





RUS Gost-R

AUS TUV CB

• GER TUV ISO13485/ISO9001/CE/CB


CHN SFDA FSC

ARG FSC EC

USA FDA & UL



YHMED 奥華 Since 1969

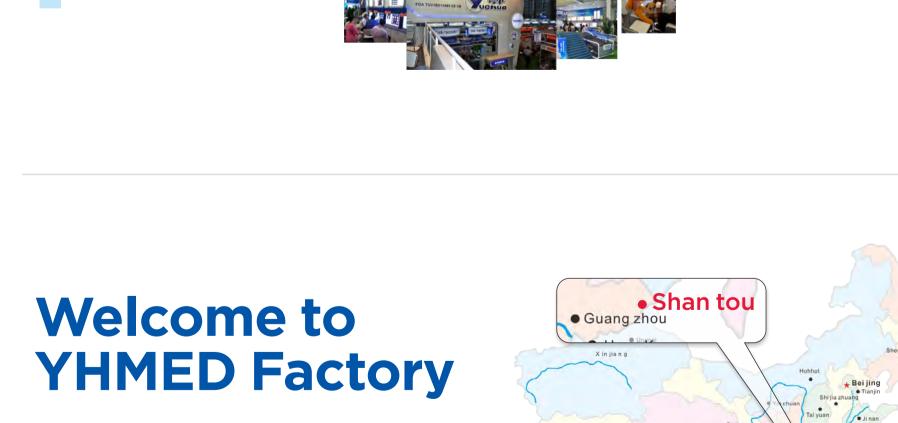

Distance

(Km/Mi)

450 / 280

350 / 220

2400 / 1500

1300 / 810

250 / 155

45 mins

X

40 mins

3 hours

2 hours

 \times

5 hours

4 hours

5 hours

4.5 hours

 \times

4 hours

3 hours

2.5 hours

 \times

1.5 hours

450 / 280 Foshan → Shantou 380 / 240 Hongkong → Shantou

Our location

Guangdong Province

Place of Departure

Guangzhou → **Shantou**

Shenzhen → Shantou

Shantou City,

Beijing → **Shantou** Shanghai → Shantou

Xiamen → Shantou

Shantou, Guangdong, China.

E-mail: export@yhmedical.com

Tel/Fax: 86-754-82525826

YHMED *** Your Best OEM & ODM Partner! Since 1969

Ms.Shirley 86-13829628128 if you have any doubts. Thank you for your kindly attention.

Please call to our Director of Overseas Marketing

@yhmed1969



Website: www.yhmedical.com yhmed.en.alibaba.com

Please check the details with our company profile video:

