

### CPAP Masks

CPAP mask with cuff. The mask cuff can be adapted to the shape of the head by the amount of air administered into the mask.

The mask has a double port. The smaller port with a 22 mm diameter is used to attach it to a tube for connection to the CPAP generator. The larger port with 30 mm is used to attach it the PEEP-valve required.


**Sterile:** Nein

Art.no.	Size
5160005521	Small/Middle
5160005520	Large

P. J. Dahlhausen & Co. GmbH  
Emil-Hoffmann-Str. 53  
50996 Köln

Tel.: (02236) 39 13 - 0  
Fax: (02236) 39 13 - 109



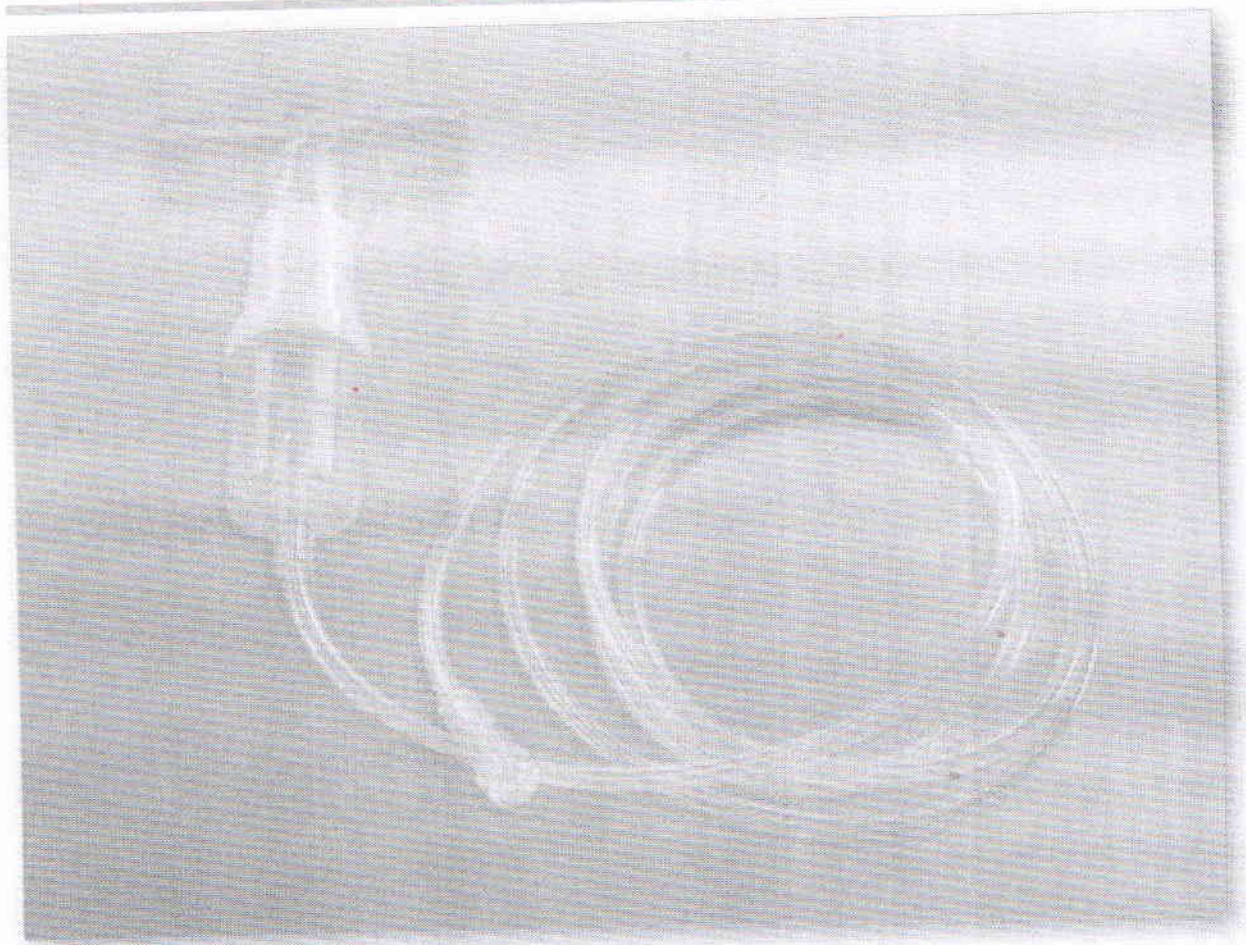

### Nebulizer Kit with standard connector

Nebulizer Kit for respiration systems is used for nebulization of medicines and their inhalation.

The nebulizer has the following outputs:

Flow rate at 3,5 bar (LPM):

4 (minimum)

**6 (recommended)**

10 (maximum)

Air pressure (bar):

1 (minimum)

**3,5 (recommended)**

6 (maximum)

MMAD (NaCl 9 g/l) ( $\mu\text{m}$ ):

2,53 (at 4 LPM)

**1,18 (at 6 LPM)**

0,88 (at 10 LPM)

Nebulized quantity (NaCl 9g/l) (mg/min)

100 (at 4 LPM)

**200 (at 6 LPM)**

340 (at 10 LPM)

Volumes of solution 2-6 ml.

Length (oxygen tube): 210 cm

Set content: Non-return valve NRV), nebulizer chamber, safety connecting tube

with two standard connectors

Sterile:

No

Art.no.	Application area
0100008905	Adult

P. J. Dahlhausen & Co. GmbH  
Emil-Hoffmann-Str. 53  
50996 Köln

Tel.: (02236) 39 13 - 0  
Fax: (02236) 39 13 - 109



DAHLHAUSEN