

RELIABLE HEARING IN ANY SITUATION

M-CORE HEARING AIDS



REXTON


OUR LATEST INNOVATIONS FOR RELIABLE HEARING



MotionCore technology



M-Core hearing aids

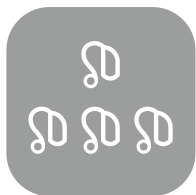

Hearing aid model overview



M-CORE B

Rechargeable BTE power, now with MotionCore

7



M-CORE R

Updated design, with Rechargeable, CROS, and Telecoil options

11



M-CORE SR

Award-winning design – now with exchangeable receivers

12



M-CORE iX

Our smallest instant-fit CIC, now with MotionCore

13

Model overview

Proven technology

14

MOTIONCORE TECHNOLOGY



Today's hearing aids analyze hearing situations solely on acoustic criteria, then classify them into a small range of predefined "acoustic environments". But when people rely on their senses to keep them safe and sound, hearing clearly in every situation is a must – especially when they're on the job and moving through the day, where environments can change quickly.

Our latest range of hearing aids use proven MotionCore technology to automatically recognize and adapt to more hearing situations than ever. This ensures users automatically hear speech clearer and get the most reliable impression of their environment. M-Core hearing aids are designed to be durable, practical and easy-to-use, so your customers can rely on Rexton in any situation.

UP TO 20 TIMES MORE HEARING SITUATIONS THAN BEFORE, AUTOMATICALLY



High Resolution Soundpro

Using a combination of inputs including real-time motion data, signal to noise ratio, and sound directionality, High Resolution Soundpro can continuously adjust the sound settings to provide the most reliable impression of the user's environment, automatically.



Motion Sensor

MotionCore technology features an integrated real-time motion sensor that detects a user's movement, and uses the data to even more accurately determine their sound situation. With this, M-Core hearing aids can automatically detect up to 127 different sound situations.



My Voice

The hearing aid recognizes the user's specific voice pattern, and integrates it naturally into what they hear. In turn, this provides the user with a natural-sounding voice that leads to high acceptance rates. And for that the user will communicate more and be more socially active, greatly improving their quality of life.



Streaming

Audio signals from smartphones, TVs or devices can be streamed wirelessly to the hearing aids; directly or via transmitter. Direct Audio Streaming is available from iOS smartphones, and Android smartphones can stream with the addition of the Smart Mic accessory.



REXTON APP: A HEARING AID WEARER'S ALL-IN-ONE SMARTPHONE COMPANION

The Rexton App allows wearers to control their hearing aids conveniently and discreetly, e.g. change programs, volume, tone.

Download the Rexton App now!



RECHARGEABLE BTE POWER

NOW WITH MOTIONCORE

Our latest generation of BTE hearing aids comes completely rechargeable, even up to our High Power models with an unbelievable runtime of up to 61 hours on a single charge. They also feature MotionCore technology, which automatically recognizes and adapts to more hearing situations than ever. With new performance levels and our first rechargeable Superpower BTE, more users can automatically get the most reliable impression of their environment than ever. M-Core B-Li hearing aids are designed to be durable, practical and easy-to-use, so your customers can rely on Rexton in any situation.

Lithium-ion power for every need

- Up to 61 hours wear time on a single full charge
- 30-minute fast charge delivers up to six hours of wear time
- Automatic drying function integrated in the optional Dry&Clean charging station



Available colors



Silver



Pearl White



Beige



Rose Gold



Gray



Dark Champagne



Black



Deep Brown



Sandy Brown



Granite



With a compact size that is comparable to a 312 battery device, the new M-CORE B-Li M combines a new, modern design with rechargeability and Bluetooth connectivity.

60/133 dB FOG/OSPL (Earhook)

60/125 dB FOG/OSPL (Life Tube)



M-CORE B-Li M

Performance Level	80, 60, 40, 30, 20
Battery	Lithium-ion battery
FOG in dB	60
Accessories	Connex Basic Charging Station Connexx Charging Station Smart Mic Smart Transmitter 2,4 New Rexton App New Life Tube and tips Color Conversion Kits
Features	Integrated motion sensor Hi-Res Soundpro Up to 120 hearing situations Direct Audio Streaming Programming via: NoahLink Wireless

Chargeable via:



Connex Basic Charging Station



Connex Charging Station

RECHARGEABLE BTES

WITH THE POWER OF MOTIONCORE

With a new smaller size, and rechargeable battery, the new M-CORE B-Li P remains a powerful device that now also offers the latest MotionCore technology.

77/135 dB FOG/OSPL (Earhook)

66/130 dBFOG/OSPL (Life Tube)

M-CORE B-LI P



Performance Level	80, 60, 40, 30, 20
Battery	Lithium-ion battery
FOG in dB	77
Accessories	Standard or Dry&Clean Charger
	Smart Mic
	Smart Transmitter 2,4
	New Rexton App
	New Life Tube and tips
	Color Conversion Kits
Features	Integrated motion sensor
	HiRes SoundPro
	Up to 120 hearing situations
	Direct Audio Streaming
	Programming via: NoahLink Wireless

Chargeable via:



B-Li P Dry&Clean Charger



Standard P Charger



Our most powerful device, the M-CORE B-Li HP is the smart choice for those with severe to profound hearing loss, combining the latest MotionCore technology with rechargeability.

82/140 dB FOG/OSPL (Earhook)



M-CORE B-LI HP

Performance Level	80, 60, 40, 30, 20
Battery	Lithium-ion battery
FOG in dB	82
Accessories	Standard or Dry&Clean Charger
	Smart Mic
	Smart Transmitter 2,4
	New Rexton App
	New Life Tube and tips
	Color Conversion Kits
Features	Integrated motion sensor (TL 8 - 16)
	HiRes SoundPro
	Up to 120 hearing situations
	Direct Audio Streaming
	Programming via: NoahLink Wireless

Chargeable via:



B-Li HP Dry&Clean Charger



Standard HP Charger

OUR MOST TRUSTED RIC

NOW WITH MOTIONCORE TECHNOLOGY

NOW IN
PERFORMANCE
LEVELS
30 & 20

Our most trusted RIC hearing aid now comes with proven MotionCore technology that automatically recognizes and adapts to more hearing situations than ever. This ensures users automatically hear speech clearer and get the most reliable impression of their environment. M-Core RIC hearing aids are designed to be durable, practical and easy-to-use – with either 312 battery or rechargeable option.

Available colors



Silver



Pearl White



Beige



Rose Gold



Gray



Dark Champagne



Black



Deep Brown



Sandy Brown



Granite



M-CORE R-LI



Performance Level	80, 60, 40 New 30, 20
Battery	Lithium-ion battery
Accessories	Charging Station
	Smart Mic
	Smart Transmitter 2,4
	Rexton App
Features	Integrated Motion Sensor
	Programmable rocker switch
	New Life Tube and tips
	Color Conversion Kits
	Programming via: NoahLink Wireless

- M-Core R-Li features a rechargeable lithium-ion battery format. The rechargeable batteries have a running time of up to 21 h incl. 5 h streaming on one full charge, and the quick charge feature delivers up to 6 hours extra running time in just 30 minutes charging.



Connexx Charging Station

M-CORE R312



Performance Level	80, 60, 40 New 30, 20
Battery	312
Accessories	Smart Mic
	Smart Transmitter 2,4
	Rexton App
Features	Integrated Motion Sensor
	Programmable rocker switch
	New Life Tube and tips
	Color Conversion Kits
	Color Conversion Kits with design integrated telephone coil (M-Core R312T)
	Programming via: NoahLink Wireless

- M-Core R312 is optionally available with an integrated Telecoil for optimal hearing where loop systems are available.



M-Core R-312T with optional telecoil

AWARD-WINNING DESIGN

NOW WITH EXCHANGEABLE RECEIVERS

NOW IN
PERFORMANCE
LEVELS
30 & 20

Our award-winning Slim-RIC design now features proven MotionCore technology that automatically recognizes and adapts to up to 20 times more hearing situations than before. This ensures users hear speech clearer and get the most reliable impression of their environment, automatically. And now with exchangeable receivers, you can easily fit M-Core SR hearing aids to more users than ever.

Choose from 5 distinct color combinations



M-CORE SR

Performance Level	80, 60, 40 New 30, 20
Battery	Lithium-ion battery
Accessories	Connexx portable charging case
	Smart Mic
	Smart Transmitter 2,4
	Rexton App
Features	Integrated Motion Sensor
	127 sound situations
	Exchangeable Receivers
	Direct Streaming
	Programming via: NoahLink Wireless

Updated portable charger for power on-the-go

M-Core SR hearing aids feature rechargeable lithium-ion batteries for hassle-free power. The included charging case enables on-the-go top-ups, for up to 4 days of use – so your customers can enjoy a relaxing weekend without having to be near a power socket.





INSTANT FIT

WITH THE HIDDEN POWER OF MOTIONCORE



Our virtually invisible CIC design now features proven MotionCore technology that automatically recognizes and adapts to up to 10 times more hearing situations than before. This ensures users hear speech clearer and get the most reliable impression of their environment, automatically. M-Core iX hearing aids are designed to be discreet, practical and easy-to-use.

M-CORE iX



Performance Level	80, 60, 40 New 30, 20
Battery	10
Accessories	Smart Key Rexton App
Features	HD Soundpro 70 sound situations Instant-Fit Sleeves Programming via: Flex connector HiPro, Connexx Air

Accessories



Smart Transmitter 2.4
Allows for direct streaming of Dolby Digital® stereo sound.











Smart Mic
Remote microphone that allows for streaming from almost any Blue-tooth-enabled phone. May also be used as a remote control.



Smart Key
Discreet remote control to change volume and programs.

A FULL RANGE OF RELIABLE HEARING OPTIONS FOR YOUR CUSTOMERS

	Custom				Instant Fit	Slim-RIC	RIC	
								
	Sterling IIC	Sterling CIC	Sterling ITC	Sterling ITE	M-CORE iX-CIC (inoX)	M-CORE SR	Emerald XS	M-CORE R312 (Emerald S)
Battery size	10	10	312	312	10/312	Li-Ion	10	312
PL80								
PL60								
PL40								
PL30	New	New	New	New	New	New		New
PL20	New	New	New	New	New	New		New
PL10								

BTE



**M-CORE
R-Li**
(Stellar RIC)

**Emerald
M**



**Mosaic
S**



**Mosaic
M**



**M-CORE
B-Li M**
(Stellar M)



**Mosaic
P**



**M-CORE
B-Li P**



**Mosaic
HP**



**M-CORE
B-Li HP**

Li-Ion	13	312	13	Li-Ion	13	Li-Ion	675	Li-Ion
● ●	● ●		● ●	● New	● ●	New	●	New
● ●	● ●		● ●	● New	● ●	New	●	New
● ●	● ●	●	● ●	● New	● ●	New	●	New
● New	●	●	● ●	● New	● ●	New	●	New
● New	●	●	● ●	● New	● ●	New	●	New
		●	●		●		●	





RELY ON REXTON

It's a tough world out there, and even tougher when your hearing starts to go. At Rexton, we won't let hearing loss hold you back. Good hearing is critical: From getting the job done, to getting home safely and being there for the people who count on you, we know what's at stake. That's why we work hard to deliver proven hearing technology in the most reliable way. We understand what reliability means in real life, and make practical, easy-to-use products you can count on. Our hearing aids have been getting the job done since 1955, so whatever life throws at you, you can rely on Rexton.

www.rexton.com

The Bluetooth® word mark and logos are owned by the Bluetooth SIG, Inc. and any use of such marks by Sivantos GmbH is under license. Other trademarks and trade names are those of their respective owners. Apple and the Apple logo are trademarks of Apple Inc., registered in the U.S. and other countries. App Store is a service mark of Apple Inc. Google Play is a trademark of Google Inc. "Made for iPhone" and "Made for iPad" mean that an electronic accessory has been designed to connect specifically to iPhone or iPad, respectively, and has been certified by the developer to meet Apple performance standards. Apple is not responsible for the operation of this device or its compliance with safety and regulatory standards. Please note that the use of this accessory with iPhone or iPad may affect wireless performance. iPad, iPhone are trademarks of Apple Inc., registered in the U.S. and other countries.

Order No. 04537-99C1-7600 © 2020 Sivantos GmbH