

Specifications 100~240VAC, 50/60Hz, 200VA Power Supply 0.1~10ml/s Flow Rate Pressure Limit 325PSI Auto Filling Rate 8ml/s Air Purge Rate 1.0~9.9ml/s Piston Rod 9.9ml/s Speed Phases 2000 Injection Programs Injection Records 2000 A: 1~200ml Volume B: 1~200ml Multiple Language User Interface

High Quality Disposables



 200ml Syringe for Contrast 200ml Syringe for Saline

antmed

Shenzhen Antmed Co., Ltd.

18 Jinhui Ave., Pingshan New District, Shenzhen 518122, China Tel: +86 755 8606 0992 E-mail: info@antmed.com

Fax: +86 755 8606 0991 http://www.antmed.com





CT Dual Injector

Leading Contrast Media Injectors Forward



CT Dual Injector

Introducing an all-new ImaStar CT dual injector, offering an extended list of product features, reliable performance, and a safe patient experience.

antmed"

Heater

Ensures contrast

temperature

maintains at ideal



Independent Dual Motor Drive System

- Dual motor drive capability for time saving loading and filling of syringes
- · Allows independent injection of contrast and saline

Dual Head Flexibility

- Simultaneous injection protocol designed for advanced cardiac CTA and Spectral CT imaging exams
- Optimization of contrast usage
- Post injection saline flush separate from contrast injection to reduce contrast wastage

Simultaneous Flow Feature

- · Allow simultaneously filling, air purging, contrast and saline injection
- · Allows the control of contrast density in both sides of the ventricles
- · Significantly reduces CT artifact
- Able to observe the right coronary arteries,
 right atrium and right ventricle in the same image

Color Touch-Screen

- 12.1 and 7" color touch-screens for easy operation
- Clear and concise information display to ensure user operational confidence

Injection Record Storage

· Stores up to 2000 injection history records



Auto Detect


· Syringe detection feature automates air purging and piston rod operation

Connect to Other Imaging Equipment

- Connect with other imaging equipment for enhanced imaging, use contrast more efficiently and minimize radiation
- Option


Real-time Pressure Curve

Display provide real-time pressure monitoring curve

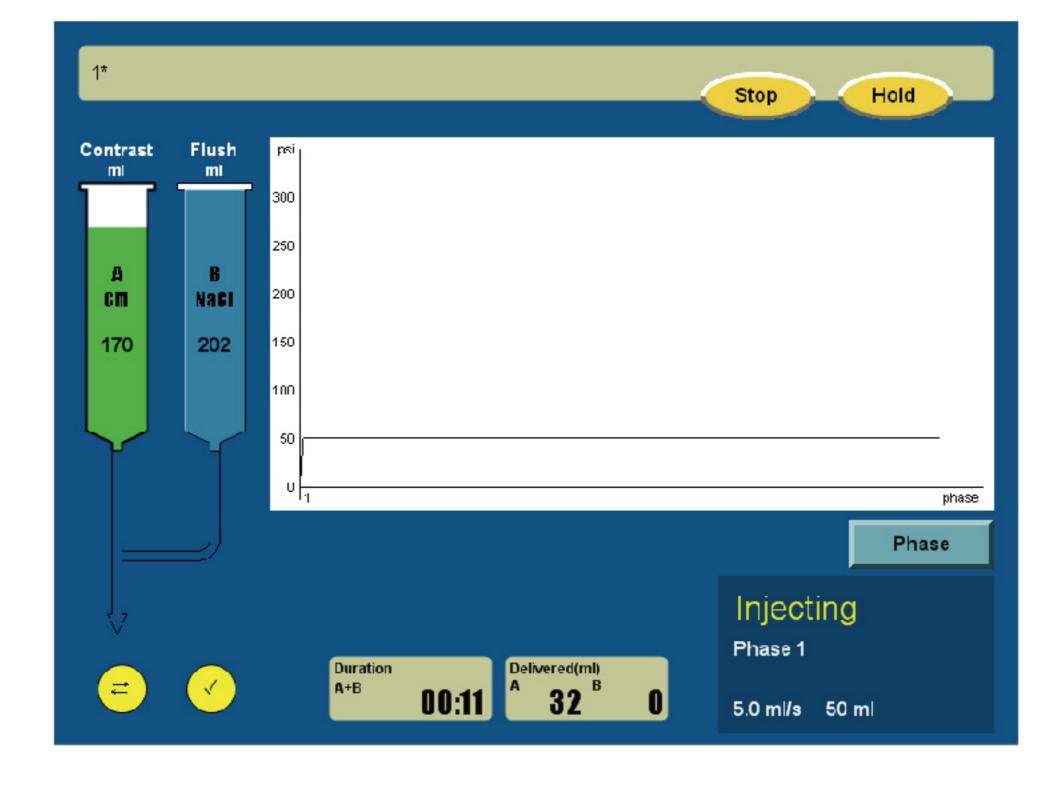

A Variety of Configuration Options

- Stores up to 2000 customized injection programs
- Program storage capacity can be increased based on user needs

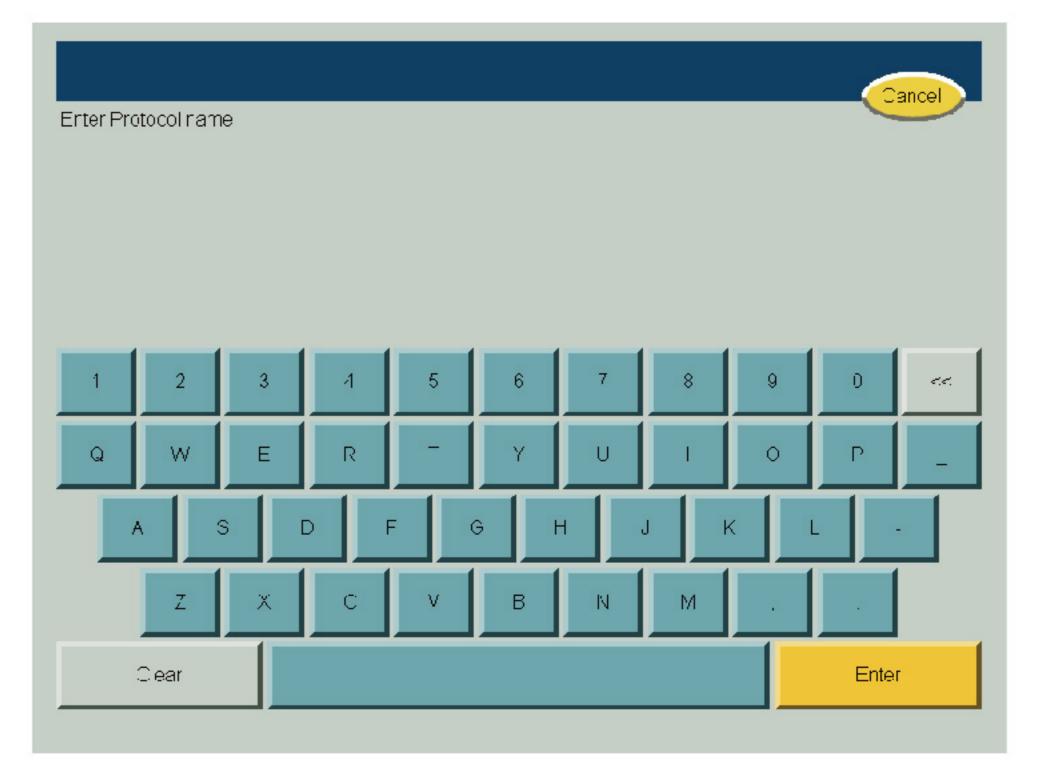

Safety

- · To ensure the safety of the injection, there is a double confirmation feature with separate prompts to confirm air purging and injection ready status before any injection
- · Automatically slows down the injection when approaching to pressure limit
- · Injection stops when pressure limit is reached. Alerts with both light and sound
- · Waterproof injector head designed to minimize damage from contrast and saline leakage
- · Ability to pause and resume while injection is underway