


Sales only to commercial customers. No sales to consumers within the meaning of §13 BGB.

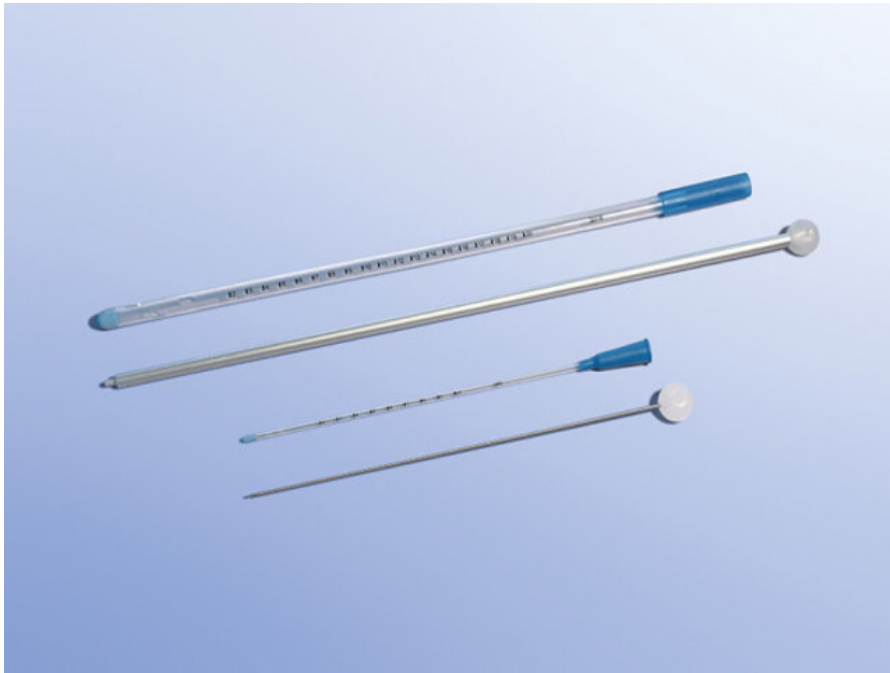

Zoom in

Trocar Catheter

The trocar catheter is inserted into the pleural gap after skin incision and puncture of the pleural cavity. The light metal stylet is removed and the remaining catheter is connected to the connecting tube/stepped connector of a thoracic drainage system.

The sizes CH 10, 12 and 16 have a funnel-shaped connector.

The sizes CH 20, 24, 28 and 32 have a plain connector.

X-ray contrastable: X-ray contrast thread

Sterile: Yes

[read more](#)

Article no.	Outer diameter	Length	Number (eyes)	minimal quantity
0200030008	8 Fg	24 cm	2	10
0200030010	10 Fg	24 cm	2	10
0200030012	12 Fg	24 cm	2	10
0200030016	16 Fg	26 cm	3	10

Article no.	Outer diameter	Length	Number (eyes)	minimal quantity
0200030020	20 Fg	40 cm	3	10
0200030024	24 Fg	40 cm	3	10
0200030028	28 Fg	40 cm	3	10
0200030032	32 Fg	40 cm	3	10