	GMD Group Medikal San. Ve Tic. A.Ş		Ref No.	DC-06
	Department	Device master file	Rev No.	03
	Damacryl and PGA	PGA (Absorbable)	Issue Date	21.07.2023

1. Product Information

1.1. General Product Name:

Sterile Single Use Absorbable PGA Surgical Suture

Product Name:

Damacryl and PGA

1.2. Code:

Suture, polyglycolic Acid: A sterile thread that is used to join the edges of a wound or incision by stitching.

GMDN CODE: 13908- Suture, polyglycolic acid

It is a sterile thread made of polyglycolic acid, which is used to stitch together the edges of a wound or incision.

CATEGORIES

03-Dental devices 10-Disposable devices


1.3. Product Images



1.4. Intended use of the product

Damacryl and PGA Surgical sutures are sterile threads that are attached to needle and are used for surgical procedures. Damacryl and PGA is indicated only for use in soft tissue approximation where wound support is required including ophthalmic surgery. Damacryl and PGA is not indicated for use in cardiovascular and neurological tissues.

The product is available in different variants. The product is available in different variants. The diameter of PGA strands varies from 0.070 mm to 0.599 mm. The length of PGA strands from 45 cm to 180 cm. The length of needles varies from 6 mm to 65 mm. The shape of needles varies from 3/8 to 1/2 circle. The different variants are listed in the PRODUCT LIST of DAMACRYL AND PGA.

	GMD Group Medikal San. Ve Tic. A.Ş		Ref No.	DC-06
	Department	Device master file	Rev No.	03
	Damacryl and PGA	PGA (Absorbable)	Issue Date	21.07.2023

Damacryl and PGA suture elicits a minimal acute inflammatory reaction in tissues, which is followed by gradual encapsulation of the suture by fibrous connective tissue. Progressive loss of tensile strength and eventual absorption of Damacryl and PGA sterile, synthetic, absorbable sutures occurs by means of hydrolysis, where the copolymer degrades to its monomers (glycolic and lactic acids), which are subsequently absorbed and metabolized in the body. Absorption begins as a loss of tensile strength without appreciable loss of mass.

The material holds approximately 65≤% of its tensile strength by two week


Surgical suture is a medical device used to hold body tissues together after an injury or surgery.

It is used for closing and covering suturable wounds under sterile conditions by initially suturing and stapling the wound.

Damacryl and PGA suture is coloured using Beige(undyed) colored multifilament braided, multifibers according to the size of the strands that meets all the requirements of the European Pharmacopoeia (EP) for Synthetic Absorbable Multifilament Sutures and the requirements of United States Pharmacopoeia (USP) Absorbable Sutures except for a slight oversize in diameter.


External appearance description

We have, undyed and Violet. It is packaged with aluminum pouch and plastic Can.

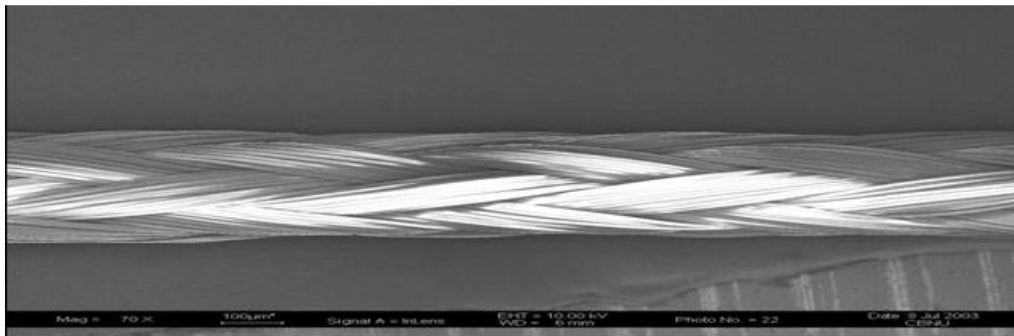

Product component

Suture is made from polyglycolic acid.

	GMD Group Medikal San. Ve Tic. A.Ş		Ref No.	DC-06
	Department	Device master file	Rev No.	03
	Damacryl and PGA	PGA (Absorbable)	Issue Date	21.07.2023



Raw Materials

Name	Raw Materials	Contents	Reference
Polymer	Polyglycolic acid	100%	-
Dye	D & C Violet No.2	under 1%	
Coating Materials	Calcium stearate	under 5%	
	Poly(caprolactone)		
Subsidiary Materials	Al pouch	Aluminium	
	PP & PE	PP & PE	
	Bobbin	-	


Product component and package size in Box

No	Contents	Components	Size	Quantity
1	MEPFIL [®] contain	PGA (Polyglycolide) suture thread	100m~1000m	12EA
3	Aluminum pouch	Al + PE	-	12EA
4	PP & PE Can	PP & PE Can	-	12EA
5	Out box	paper	-	12EA
6	Label	paper	-	12EA

Product appearance

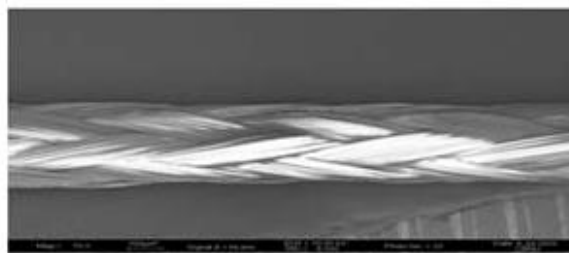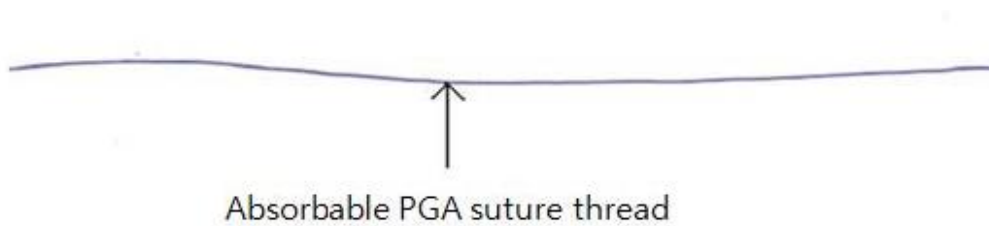


Enlarged picture (SEM)


	GMD Group Medikal San. Ve Tic. A.Ş		Ref No.	DC-06
	Department	Device master file	Rev No.	03
	Damacryl and PGA	PGA (Absorbable)	Issue Date	21.07.2023

Product codes and Description:


ORIGINAL PRODUCT BARCODE	CODE DESCRIPTION	PRODUCT REFERENCE CODE	BRAND NAME	GMDN	CLASS	SH LI (M
8681665009144	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 150cm USP SIZE 0 (VIOLET) WITHOUT NEEDLE	DC0-15	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8681677073324	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 150cm USP SIZE 2 (VIOLET) WITHOUT NEEDLE	DC02-15	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8681677073614	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 2PCS X 150cm USP SIZE 2 (VIOLET) WITHOUT NEEDLE	DC02-2X150	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8681677075656	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 2PCS X 75cm USP SIZE 2 (VIOLET) WITHOUT NEEDLE	DC02-2X75	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8681677074451	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 5PCS X 75cm USP SIZE 2 (VIOLET) WITHOUT NEEDLE	DC02-5X75	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60

	GMD Group Medikal San. Ve Tic. A.Ş		Ref No.	DC-06
	Department	Device master file	Rev No.	03
	Damacryl and PGA	PGA (Absorbable)	Issue Date	21.07.2023


8681677073621	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 7PCS X 50cm USP SIZE 2 (VIOLET) WITHOUT NEEDLE	DC02-7X50	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8681677071672	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 75cm USP SIZE 2 (VIOLET) 40mm 3/8 CIRCLE REVERSE CUTTING NEEDLE	DC02K40	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8681677074772	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 90cm USP SIZE 2 (VIOLET) 40mm 3/8 CIRCLE REVERSE CUTTING NEEDLE	DC02K40-90	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8681677072518	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 75cm USP SIZE 2 (VIOLET) 45mm 3/8 CIRCLE REVERSE CUTTING NEEDLE	DC02K45	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8681677074468	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 75cm USP SIZE 2 (VIOLET) 48mm 3/8 CIRCLE REVERSE CUTTING NEEDLE	DC02K48	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8681677073850	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 120cm USP SIZE 2 (VIOLET) 48mm 3/8 CIRCLE REVERSE CUTTING NEEDLE	DC02K48-120	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8681677075489	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 90cm USP SIZE 2 (VIOLET) 48mm 3/8 CIRCLE REVERSE CUTTING NEEDLE	DC02K48-90	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60

	GMD Group Medikal San. Ve Tic. A.Ş		Ref No.	DC-06
	Department	Device master file	Rev No.	03
	Damacryl and PGA	PGA (Absorbable)	Issue Date	21.07.2023


8681677072471	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 75cm USP SIZE 2 (VIOLET) 50mm 3/8 CIRCLE REVERSE CUTTING NEEDLE	DC02K50	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8681677075434	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 100cm USP SIZE 2 (VIOLET) 50mm 3/8 CIRCLE REVERSE CUTTING NEEDLE	DC02K50-100	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8681665009229	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 75cm USP SIZE 2 (VIOLET) 40mm 1/2 CIRCLE TAPER POINT NEEDLE	DC02M40	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8681677074130	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 120cm USP SIZE 2 (VIOLET) 40mm 1/2 CIRCLE TAPER POINT NEEDLE	DC02M40-120	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8681677073812	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 150cm USP SIZE 2 (VIOLET) 40mm 1/2 CIRCLE TAPER POINT NEEDLE	DC02M40-150	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8681677070453	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 90cm USP SIZE 2 (VIOLET) 40mm 1/2 CIRCLE TAPER POINT NEEDLE	DC02M40-90	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8681677071658	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 75cm USP SIZE 2 (VIOLET) 45mm 1/2 CIRCLE TAPER POINT NEEDLE	DC02M45	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60

	GMD Group Medikal San. Ve Tic. A.Ş		Ref No.	DC-06
	Department	Device master file	Rev No.	03
	Damacryl and PGA	PGA (Absorbable)	Issue Date	21.07.2023


8681677072501	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 75cm USP SIZE 2 (VIOLET) 48mm 1/2 CIRCLE TAPER POINT NEEDLE	DC02M48	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8681677074413	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 100cm USP SIZE 2 (VIOLET) 48mm 1/2 CIRCLE TAPER POINT NEEDLE	DC02M48-100	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8681677073829	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 120cm USP SIZE 2 (VIOLET) 48mm 1/2 CIRCLE TAPER POINT NEEDLE	DC02M48-120	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8681677072570	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 90cm USP SIZE 2 (VIOLET) 48mm 1/2 CIRCLE TAPER POINT NEEDLE	DC02M48-90	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8681677071221	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 75cm USP SIZE 2 (VIOLET) 50mm 1/2 CIRCLE TAPER POINT NEEDLE	DC02M50	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8681677073997	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 100cm USP SIZE 2 (VIOLET) 50mm 1/2 CIRCLE TAPER POINT NEEDLE	DC02M50-100	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8681677074963	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 150cm USP SIZE 2 (VIOLET) 50mm 1/2 CIRCLE TAPER POINT NEEDLE	DC02M50-150	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60

	GMD Group Medikal San. Ve Tic. A.Ş		Ref No.	DC-06
	Department	Device master file	Rev No.	03
	Damacryl and PGA	PGA (Absorbable)	Issue Date	21.07.2023


8681677071665	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 90cm USP SIZE 2 (VIOLET) 50mm 1/2 CIRCLE TAPER POINT NEEDLE	DC02M50-90	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8681677073652	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 2PCS X 150cm USP SIZE 0 (VIOLET) WITHOUT NEEDLE	DC0-2X150	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8681677075649	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 2PCS X 75cm USP SIZE 0 (VIOLET) WITHOUT NEEDLE	DC0-2X75	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8681677073669	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 7PCS X 50cm USP SIZE 0 (VIOLET) WITHOUT NEEDLE	DC0-7X50	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8681677075496	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 75cm USP SIZE 0 (VIOLET) 30mm 3/8 CIRCLE REVERSE CUTTING NEEDLE	DC0K30	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8681677073522	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 90cm USP SIZE 0 (VIOLET) 30mm 3/8 CIRCLE REVERSE CUTTING NEEDLE	DC0K30-90	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8681677070859	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 75cm USP SIZE 0 (VIOLET) 36mm 3/8 CIRCLE REVERSE CUTTING NEEDLE	DC0K36	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60

	GMD Group Medikal San. Ve Tic. A.Ş		Ref No.	DC-06
	Department	Device master file	Rev No.	03
	Damacryl and PGA	PGA (Absorbable)	Issue Date	21.07.2023


8681677074789	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 90cm USP SIZE 0 (VIOLET) 36mm 3/8 CIRCLE REVERSE CUTTING NEEDLE	DC0K36-90	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8681665009151	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 75cm USP SIZE 0 (VIOLET) 40mm 3/8 CIRCLE REVERSE CUTTING NEEDLE	DC0K40	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8681677075687	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 100cm USP SIZE 0 (VIOLET) 40mm 3/8 CIRCLE REVERSE CUTTING NEEDLE	DC0K40-100	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8681677074178	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 120cm USP SIZE 0 (VIOLET) 40mm 3/8 CIRCLE REVERSE CUTTING NEEDLE	DC0K40-120	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8681677074123	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 90cm USP SIZE 0 (VIOLET) 40mm 3/8 CIRCLE REVERSE CUTTING NEEDLE	DC0K40-90	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8681677072495	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 75cm USP SIZE 0 (VIOLET) 45mm 3/8 CIRCLE REVERSE CUTTING NEEDLE	DC0K45	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8681677074260	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 75cm USP SIZE 0 (VIOLET) 48mm 3/8 CIRCLE REVERSE CUTTING NEEDLE	DC0K48	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60

	GMD Group Medikal San. Ve Tic. A.Ş		Ref No.	DC-06
	Department	Device master file	Rev No.	03
	Damacryl and PGA	PGA (Absorbable)	Issue Date	21.07.2023


8681677072464	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 75cm USP SIZE 0 (VIOLET) 50mm 3/8 CIRCLE REVERSE CUTTING NEEDLE	DC0K50	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8681677074109	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 100cm USP SIZE 0 (VIOLET) 50mm 3/8 CIRCLE REVERSE CUTTING NEEDLE	DC0K50-100	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8681665009168	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 75cm USP SIZE 0 (VIOLET) 30mm 1/2 CIRCLE TAPER POINT NEEDLE	DC0M30	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8681677073843	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 120cm USP SIZE 0 (VIOLET) 30mm 1/2 CIRCLE TAPER POINT NEEDLE	DC0M30-120	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8681677072877	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 90cm USP SIZE 0 (VIOLET) 30mm 1/2 CIRCLE TAPER POINT NEEDLE	DC0M30-90	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8681677071603	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 75cm USP SIZE 0 (VIOLET) 35mm 1/2 CIRCLE TAPER POINT NEEDLE	DC0M35	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8681677072556	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 90cm USP SIZE 0 (VIOLET) 35mm 1/2 CIRCLE TAPER POINT NEEDLE	DC0M35-90	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60

	GMD Group Medikal San. Ve Tic. A.Ş		Ref No.	DC-06
	Department	Device master file	Rev No.	03
	Damacryl and PGA	PGA (Absorbable)	Issue Date	21.07.2023


8681677070118	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 75cm USP SIZE 0 (VIOLET) 36mm 1/2 CIRCLE TAPER POINT NEEDLE	DC0M36	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8681677074161	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 120cm USP SIZE 0 (VIOLET) 36mm 1/2 CIRCLE TAPER POINT NEEDLE	DC0M36-120	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8681677072815	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 90cm USP SIZE 0 (VIOLET) 36mm 1/2 CIRCLE TAPER POINT NEEDLE	DC0M36-90	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8681665009175	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 75cm USP SIZE 0 (VIOLET) 40mm 1/2 CIRCLE TAPER POINT NEEDLE	DC0M40	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8681677072631	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 100cm USP SIZE 0 (VIOLET) 40mm 1/2 CIRCLE TAPER POINT NEEDLE	DC0M40-100	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8681677070279	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 90cm USP SIZE 0 (VIOLET) 40mm 1/2 CIRCLE TAPER POINT NEEDLE	DC0M40-90	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8681677072426	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 75cm USP SIZE 0 (VIOLET) 45mm 1/2 CIRCLE TAPER POINT NEEDLE	DC0M45	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60

	GMD Group Medikal San. Ve Tic. A.Ş		Ref No.	DC-06
	Department	Device master file	Rev No.	03
	Damacryl and PGA	PGA (Absorbable)	Issue Date	21.07.2023


8681677074345	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 90cm USP SIZE 0 (VIOLET) 45mm 1/2 CIRCLE TAPER POINT NEEDLE	DC0M45-90	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8681677073607	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 75cm USP SIZE 0 (VIOLET) 48mm 1/2 CIRCLE TAPER POINT NEEDLE	DC0M48	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8681677072198	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 75cm USP SIZE 0 (VIOLET) 50mm 1/2 CIRCLE TAPER POINT NEEDLE	DC0M50	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8681677072655	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 100cm USP SIZE 0 (VIOLET) 50mm 1/2 CIRCLE TAPER POINT NEEDLE	DC0M50-100	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8681677071610	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 90cm USP SIZE 0 (VIOLET) 50mm 1/2 CIRCLE TAPER POINT NEEDLE	DC0M50-90	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8681677073355	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 150cm USP SIZE 1 (VIOLET) WITHOUT NEEDLE	DC1-15	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8681677073638	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 2PCS X 150cm USP SIZE 1 (VIOLET) WITHOUT NEEDLE	DC1-2X150	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60

	GMD Group Medikal San. Ve Tic. A.Ş		Ref No.	DC-06
	Department	Device master file	Rev No.	03
	Damacryl and PGA	PGA (Absorbable)	Issue Date	21.07.2023


8681677075632	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 2PCS X 75cm USP SIZE 1 (VIOLET) WITHOUT NEEDLE	DC1-2X75	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8681677073348	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 5PCS X 75cm USP SIZE 1 (VIOLET DYED) WITHOUT NEEDLE	DC1-7X5	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8681677073645	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 7PCS X 50cm USP SIZE 1 (VIOLET) WITHOUT NEEDLE	DC1-7X50	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8681677075663	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 75cm USP SIZE 1 (VIOLET) 30mm 3/8 CIRCLE REVERSE CUTTING NEEDLE	DC1K30	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8681677074338	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 75cm USP SIZE 1 (VIOLET) 35mm 3/8 CIRCLE REVERSE CUTTING NEEDLE	DC1K35	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8681677073300	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 75cm USP SIZE 1 (VIOLET) 36mm 3/8 CIRCLE REVERSE CUTTING NEEDLE	DC1K36	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8681665009182	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 75cm USP SIZE 1 (VIOLET) 37mm 3/8 CIRCLE REVERSE CUTTING NEEDLE	DC1K37	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60

	GMD Group Medikal San. Ve Tic. A.Ş		Ref No.	DC-06
	Department	Device master file	Rev No.	03
	Damacryl and PGA	PGA (Absorbable)	Issue Date	21.07.2023


8681665009199	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 75cm USP SIZE 1 (VIOLET) 40mm 3/8 CIRCLE REVERSE CUTTING NEEDLE	DC1K40	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8681677074154	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 120cm USP SIZE 1 (VIOLET) 40mm 3/8 CIRCLE REVERSE CUTTING NEEDLE	DC1K40-120	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8681677072709	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 90cm USP SIZE 1 (VIOLET) 40mm 3/8 CIRCLE REVERSE CUTTING NEEDLE	DC1K40-90	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8681677072488	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 75cm USP SIZE 1 (VIOLET) 45mm 3/8 CIRCLE REVERSE CUTTING NEEDLE	DC1K45	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8681677074253	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 75cm USP SIZE 1 (VIOLET) 48mm 3/8 CIRCLE REVERSE CUTTING NEEDLE	DC1K48	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8681677074796	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 90cm USP SIZE 1 (VIOLET) 48mm 3/8 CIRCLE REVERSE CUTTING NEEDLE	DC1K48-90	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8681677072457	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 75cm USP SIZE 1 (VIOLET) 50mm 3/8 CIRCLE REVERSE CUTTING NEEDLE	DC1K50	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60

	GMD Group Medikal San. Ve Tic. A.Ş		Ref No.	DC-06
	Department	Device master file	Rev No.	03
	Damacryl and PGA	PGA (Absorbable)	Issue Date	21.07.2023


8681677074055	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 100cm USP SIZE 1 (VIOLET) 50mm 3/8 CIRCLE REVERSE CUTTING NEEDLE	DC1K50-100	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8681677075762	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 75cm USP SIZE 1 (VIOLET) 40mm 1/2 CIRCLE NORMAL CUTTING EDGE NEEDLE	DC1KDN40	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8681665009205	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 75cm USP SIZE 1 (VIOLET) 30mm 1/2 CIRCLE TAPER POINT NEEDLE	DC1M30	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8681677074147	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 120cm USP SIZE 1 (VIOLET) 30mm 1/2 CIRCLE TAPER POINT NEEDLE	DC1M30-120	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8681677073515	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 90cm USP SIZE 1 (VIOLET) 30mm 1/2 CIRCLE TAPER POINT NEEDLE	DC1M30-90	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8681677075502	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 75cm USP SIZE 1 (VIOLET) 35mm 1/2 CIRCLE TAPER POINT NEEDLE	DC1M35	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8681677075519	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 120cm USP SIZE 1 (VIOLET) 35mm 1/2 CIRCLE TAPER POINT NEEDLE	DC1M35-120	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60

	GMD Group Medikal San. Ve Tic. A.Ş		Ref No.	DC-06
	Department	Device master file	Rev No.	03
	Damacryl and PGA	PGA (Absorbable)	Issue Date	21.07.2023


8681677073270	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 90cm USP SIZE 1 (VIOLET) 35mm 1/2 CIRCLE TAPER POINT NEEDLE	DC1M35-90	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8681665009212	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 75cm USP SIZE 1 (VIOLET) 40mm 1/2 CIRCLE TAPER POINT NEEDLE	DC1M40	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8681677072679	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 100cm USP SIZE 1 (VIOLET) 40mm 1/2 CIRCLE TAPER POINT NEEDLE	DC1M40-100	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8681677073799	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 120cm USP SIZE 1 (VIOLET) 40mm 1/2 CIRCLE TAPER POINT NEEDLE	DC1M40-120	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8681677070286	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 90cm USP SIZE 1 (VIOLET) 40mm 1/2 CIRCLE TAPER POINT NEEDLE	DC1M40-90	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8681677071627	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 75cm USP SIZE 1 (VIOLET) 45mm 1/2 CIRCLE TAPER POINT NEEDLE	DC1M45	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8681677072624	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 100cm USP SIZE 1 (VIOLET) 45mm 1/2 CIRCLE TAPER POINT NEEDLE	DC1M45-100	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60

	GMD Group Medikal San. Ve Tic. A.Ş		Ref No.	DC-06
	Department	Device master file	Rev No.	03
	Damacryl and PGA	PGA (Absorbable)	Issue Date	21.07.2023


8681677074314	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 90cm USP SIZE 1 (VIOLET) 45mm 1/2 CIRCLE TAPER POINT NEEDLE	DC1M45-90	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8681677072440	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 75cm USP SIZE 1 (VIOLET) 48mm 1/2 CIRCLE TAPER POINT NEEDLE	DC1M48	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8681677072617	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 100cm USP SIZE 1 (VIOLET) 48mm 1/2 CIRCLE TAPER POINT NEEDLE	DC1M48-100	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8681677072693	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 150cm USP SIZE 1 (VIOLET) 48mm 1/2 CIRCLE TAPER POINT DOUBLE NEEDLE	DC1M48-150-D	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8681677072563	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 90cm USP SIZE 1 (VIOLET) 48mm 1/2 CIRCLE TAPER POINT NEEDLE	DC1M48-90	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8681677071634	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 75cm USP SIZE 1 (VIOLET) 50mm 1/2 CIRCLE TAPER POINT NEEDLE	DC1M50	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8681677073294	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 100cm USP SIZE 1 (VIOLET) 50mm 1/2 CIRCLE TAPER POINT NEEDLE	DC1M50-100	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60

	GMD Group Medikal San. Ve Tic. A.Ş		Ref No.	DC-06
	Department	Device master file	Rev No.	03
	Damacryl and PGA	PGA (Absorbable)	Issue Date	21.07.2023


8681677071641	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 90cm USP SIZE 1 (VIOLET) 50mm 1/2 CIRCLE TAPER POINT NEEDLE	DC1M50-90	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8681677073805	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 150cm USP SIZE 1 (VIOLET) 40mm 1/2 CIRCLE TAPER POINT NEEDLE (LOOPED ON)	DC1ML40-150	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8681665009236	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 150cm USP SIZE 2-0 (VIOLET DYED) WITHOUT NEEDLE	DC2-15	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8681677073362	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 5PCS X 75cm USP SIZE 2-0 (VIOLET) WITHOUT NEEDLE	DC2-7X5	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8681665009243	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 75cm USP SIZE 2-0 (VIOLET) 26mm 3/8 CIRCLE REVERSE CUTTING NEEDLE	DC2K26	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8681665009304	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 75cm USP SIZE 2-0 (UNDYED) 26mm 3/8 CIRCLE REVERSE CUTTING NEEDLE	DC2K26U	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8681665009250	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 75cm USP SIZE 2-0 (VIOLET) 30mm 3/8 CIRCLE REVERSE CUTTING NEEDLE	DC2K30	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60

	GMD Group Medikal San. Ve Tic. A.Ş		Ref No.	DC-06
	Department	Device master file	Rev No.	03
	Damacryl and PGA	PGA (Absorbable)	Issue Date	21.07.2023


8681677073836	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 100cm USP SIZE 2-0 (VIOLET) 30mm 3/8 CIRCLE REVERSE CUTTING NEEDLE	DC2K30-100	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8681677073539	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 90cm USP SIZE 2-0 (VIOLET) 30mm 3/8 CIRCLE REVERSE CUTTING NEEDLE	DC2K30-90	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8681677073744	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 75cm USP SIZE 2-0 (VIOLET) 35mm 3/8 CIRCLE REVERSE CUTTING NEEDLE	DC2K35	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8681677074208	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 120cm USP SIZE 2-0 (VIOLET) 36mm 3/8 CIRCLE REVERSE CUTTING NEEDLE	DC2K36-120	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8681677074215	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 90cm USP SIZE 2-0 (VIOLET) 36mm 3/8 CIRCLE REVERSE CUTTING NEEDLE	DC2K36-90	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8681677074680	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 75cm USP SIZE 2-0 (VIOLET) 40mm 3/8 CIRCLE REVERSE CUTTING NEEDLE	DC2K40	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8681677071597	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 75cm USP SIZE 2-0 (VIOLET) 19mm 1/2 CIRCLE TAPER POINT NEEDLE	DC2M19	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60

	GMD Group Medikal San. Ve Tic. A.Ş		Ref No.	DC-06
	Department	Device master file	Rev No.	03
	Damacryl and PGA	PGA (Absorbable)	Issue Date	21.07.2023


8681677071566	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 75cm USP SIZE 2-0 (VIOLET) 20mm 1/2 CIRCLE TAPER POINT NEEDLE	DC2M20	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8681665009267	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 75cm USP SIZE 2-0 (VIOLET) 26mm 1/2 CIRCLE TAPER POINT NEEDLE	DC2M26	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8681677073775	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 100cm USP SIZE 2-0 (VIOLET) 26mm 1/2 CIRCLE TAPER POINT NEEDLE	DC2M26-100	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8681677074192	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 120cm USP SIZE 2-0 (VIOLET) 26mm 1/2 CIRCLE TAPER POINT NEEDLE	DC2M26-120	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8681665009274	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 75cm USP SIZE 2-0 (VIOLET) 30mm 1/2 CIRCLE TAPER POINT NEEDLE	DC2M30	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8681677073782	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 100cm USP SIZE 2-0 (VIOLET) 30mm 1/2 CIRCLE TAPER POINT NEEDLE	DC2M30-100	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8681677074185	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 120cm USP SIZE 2-0 (VIOLET) 30mm 1/2 CIRCLE TAPER POINT NEEDLE	DC2M30-120	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60

	GMD Group Medikal San. Ve Tic. A.Ş		Ref No.	DC-06
	Department	Device master file	Rev No.	03
	Damacryl and PGA	PGA (Absorbable)	Issue Date	21.07.2023


8681677072884	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 90cm USP SIZE 2-0 (VIOLET) 30mm 1/2 CIRCLE TAPER POINT NEEDLE	DC2M30-90	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8681677071573	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 75cm USP SIZE 2-0 (VIOLET) 35mm 1/2 CIRCLE TAPER POINT NEEDLE	DC2M35	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8681677072907	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 90cm USP SIZE 2-0 (VIOLET) 35mm 1/2 CIRCLE TAPER POINT NEEDLE	DC2M35-90	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8681665009281	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 75cm USP SIZE 2-0 (VIOLET) 36mm 1/2 CIRCLE TAPER POINT NEEDLE	DC2M36	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8681677074765	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 90cm USP SIZE 2-0 (VIOLET) 36mm 1/2 CIRCLE TAPER POINT NEEDLE	DC2M36-90	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8681665009298	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 75cm USP SIZE 2-0 (VIOLET) 40mm 1/2 CIRCLE TAPER POINT NEEDLE	DC2M40	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8681677074000	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 100cm USP SIZE 2-0 (VIOLET) 40mm 1/2 CIRCLE TAPER POINT NEEDLE	DC2M40-100	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60

	GMD Group Medikal San. Ve Tic. A.Ş		Ref No.	DC-06
	Department	Device master file	Rev No.	03
	Damacryl and PGA	PGA (Absorbable)	Issue Date	21.07.2023


8681677072549	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 90cm USP SIZE 2-0 (VIOLET) 40mm 1/2 CIRCLE TAPER POINT NEEDLE	DC2M40-90	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8681677074437	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 75cm USP SIZE 2-0 (VIOLET) 48mm 1/2 CIRCLE TAPER POINT NEEDLE	DC2M48	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8681677071580	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 75cm USP SIZE 2-0 (VIOLET) 50mm 1/2 CIRCLE TAPER POINT NEEDLE	DC2M50	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8681677072433	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 75cm USP SIZE 2-0 (VIOLET) 60mm STRAIGHT REVERSE CUTTING NEEDLE	DC2SK60	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8681677070132	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 2PCS X 75cm USP SIZE 3-0 (VIOLET) WITHOUT NEEDLE	DC3-7X2	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8681677071320	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 6PCS X 75cm USP SIZE 3-0 (VIOLET) WITHOUT NEEDLE	DC3-7X6	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8681677074239	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 120cm USP SIZE 3-0 (VIOLET) 16mm 3/8 CIRCLE REVERSE CUTTING NEEDLE	DC3K16-120	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60

	GMD Group Medikal San. Ve Tic. A.Ş		Ref No.	DC-06
	Department	Device master file	Rev No.	03
	Damacryl and PGA	PGA (Absorbable)	Issue Date	21.07.2023


8681677071559	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 75cm USP SIZE 3-0 (VIOLET) 19mm 3/8 CIRCLE REVERSE CUTTING NEEDLE	DC3K19	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8681677071368	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 75cm USP SIZE 3-0 (VIOLET) 22mm 3/8 CIRCLE REVERSE CUTTING NEEDLE	DC3K22	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8681665009342	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 75cm USP SIZE 3-0 (VIOLET) 26mm 3/8 CIRCLE REVERSE CUTTING NEEDLE	DC3K26	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8681665009380	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 75cm USP SIZE 3-0 (UNDYED) 26mm 3/8 CIRCLE REVERSE CUTTING NEEDLE	DC3K26U	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8681677075618	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 75cm USP SIZE 3-0 (VIOLET) 30mm 3/8 CIRCLE REVERSE CUTTING NEEDLE	DC3K30	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8681665009335	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 75cm USP SIZE 3-0 (VIOLET) 22mm 3/8 CIRCLE NORMAL CUTTING NEEDLE	DC3KD22	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8681665009311	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 45cm USP SIZE 3-0 (VIOLET) 22mm 3/8 CIRCLE NORMAL CUTTING NEEDLE	DC3KD22-4	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60

	GMD Group Medikal San. Ve Tic. A.Ş		Ref No.	DC-06
	Department	Device master file	Rev No.	03
	Damacryl and PGA	PGA (Absorbable)	Issue Date	21.07.2023


8681677072310	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 75cm USP SIZE 3-0 (VIOLET) 18mm 1/2 CIRCLE TAPER POINT NEEDLE	DC3M18	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8681665009359	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 75cm USP SIZE 3-0 (VIOLET) 20mm 1/2 CIRCLE TAPER POINT NEEDLE	DC3M20	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8681677073751	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 100cm USP SIZE 3-0 (VIOLET) 20mm 1/2 CIRCLE TAPER POINT NEEDLE	DC3M20-100	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8681665009366	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 75cm USP SIZE 3-0 (VIOLET) 22mm 1/2 CIRCLE TAPER POINT NEEDLE	DC3M22	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8681677071535	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 75cm USP SIZE 3-0 (VIOLET) 25mm 1/2 CIRCLE TAPER POINT NEEDLE	DC3M25	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8681665009373	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 75cm USP SIZE 3-0 (VIOLET) 26mm 1/2 CIRCLE TAPER POINT NEEDLE	DC3M26	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8681677073768	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 100cm USP SIZE 3-0 (VIOLET) 26mm 1/2 CIRCLE TAPER POINT NEEDLE	DC3M26-100	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60

	GMD Group Medikal San. Ve Tic. A.Ş		Ref No.	DC-06
	Department	Device master file	Rev No.	03
	Damacryl and PGA	PGA (Absorbable)	Issue Date	21.07.2023


8681677074222	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 120cm USP SIZE 3-0 (VIOLET) 26mm 1/2 CIRCLE TAPER POINT NEEDLE	DC3M26-120	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8681677074956	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 90cm USP SIZE 3-0 (VIOLET) 26mm 1/2 CIRCLE TAPER POINT NEEDLE	DC3M26-90	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8681677071542	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 75cm USP SIZE 3-0 (VIOLET) 30mm 1/2 CIRCLE TAPER POINT NEEDLE	DC3M30	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8681677074017	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 100cm USP SIZE 3-0 (VIOLET) 30mm 1/2 CIRCLE TAPER POINT NEEDLE	DC3M30-100	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8681677072938	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 75cm USP SIZE 3-0 (VIOLET) 60mm STRAIGHT REVERSE CUTTING NEEDLE	DC3SK60	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8681677074369	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 150cm USP SIZE 4-0 (VIOLET) WITHOUT NEEDLE	DC4-150	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8681677070545	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 75cm USP SIZE 4-0 (VIOLET) 16mm 3/8 CIRCLE REVERSE CUTTING NEEDLE	DC4K16	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60

	GMD Group Medikal San. Ve Tic. A.Ş		Ref No.	DC-06
	Department	Device master file	Rev No.	03
	Damacryl and PGA	PGA (Absorbable)	Issue Date	21.07.2023


8681677071528	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 75cm USP SIZE 4-0 (VIOLET) 19mm 3/8 CIRCLE REVERSE CUTTING NEEDLE	DC4K19	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8681677072945	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 75cm USP SIZE 4-0 (VIOLET) 20mm 3/8 CIRCLE REVERSE CUTTING NEEDLE	DC4K20	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8681677073232	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 75cm USP SIZE 4-0 (VIOLET) 26mm 3/8 CIRCLE REVERSE CUTTING NEEDLE	DC4K26	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8681665009410	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 75cm USP SIZE 4-0 (UNDYED) 26mm 3/8 CIRCLE REVERSE CUTTING NEEDLE	DC4K26U	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8681665009397	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 75cm USP SIZE 4-0 (VIOLET) 16mm 3/8 CIRCLE NORMAL CUTTING NEEDLE	DC4KD16	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8681677073546	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 75cm USP SIZE 4-0 (VIOLET) 12mm 1/2 CIRCLE TAPER POINT NEEDLE	DC4M12	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8681677070811	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 75cm USP SIZE 4-0 (VIOLET) 16mm 1/2 CIRCLE TAPER POINT NEEDLE	DC4M16	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60

	GMD Group Medikal San. Ve Tic. A.Ş		Ref No.	DC-06
	Department	Device master file	Rev No.	03
	Damacryl and PGA	PGA (Absorbable)	Issue Date	21.07.2023


8681677074246	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 120cm USP SIZE 4-0 (VIOLET) 16mm 1/2 CIRCLE TAPER POINT NEEDLE	DC4M16-120	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8681677073867	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 75cm USP SIZE 4-0 (VIOLET) 17mm 1/2 CIRCLE TAPER POINT NEEDLE	DC4M17	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8681665009403	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 75cm USP SIZE 4-0 (VIOLET) 20mm 1/2 CIRCLE TAPER POINT NEEDLE	DC4M20	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8681677070538	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 75cm USP SIZE 4-0 (VIOLET) 22mm 1/2 CIRCLE TAPER POINT NEEDLE	DC4M22	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8681677073317	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 75cm USP SIZE 4-0 (VIOLET) 26mm 1/2 CIRCLE TAPER POINT NEEDLE	DC4M26	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8681677074376	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 150cm USP SIZE 5-0 (VIOLET) WITHOUT NEEDLE	DC5-150	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8681677071467	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 75cm USP SIZE 5-0 (VIOLET) 13mm 3/8 CIRCLE REVERSE CUTTING NEEDLE	DC5K13	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60

	GMD Group Medikal San. Ve Tic. A.Ş		Ref No.	DC-06
	Department	Device master file	Rev No.	03
	Damacryl and PGA	PGA (Absorbable)	Issue Date	21.07.2023


8681665009427	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 75cm USP SIZE 5-0 (VIOLET) 16mm 3/8 CIRCLE REVERSE CUTTING NEEDLE	DC5K16	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8681677074802	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 75cm USP SIZE 5-0 (VIOLET) 17mm 3/8 CIRCLE REVERSE CUTTING NEEDLE	DC5K17	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8681677071511	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 75cm USP SIZE 5-0 (VIOLET) 13mm 1/2 CIRCLE TAPER POINT NEEDLE	DC5M13	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8681665009434	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 75cm USP SIZE 5-0 (VIOLET) 16mm 1/2 CIRCLE TAPER POINT NEEDLE	DC5M16	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8681677074383	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 150cm USP SIZE 6-0 (VIOLET) WITHOUT NEEDLE	DC6-150	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8681677070569	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 75cm USP SIZE 6-0 (VIOLET) 13mm 3/8 CIRCLE REVERSE CUTTING NEEDLE	DC6K13	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8681677070774	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 75cm USP SIZE 6-0 (VIOLET) 13mm 1/2 CIRCLE TAPER POINT NEEDLE	DC6M13	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60

	GMD Group Medikal San. Ve Tic. A.Ş		Ref No.	DC-06
	Department	Device master file	Rev No.	03
	Damacryl and PGA	PGA (Absorbable)	Issue Date	21.07.2023


8681677074482	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 75cm USP SIZE 6-0 (VIOLET) 11mm 1/2 CIRCLE TAPER POINT NEEDLE	DC6M11	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8681677075854	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 120cm USP SIZE 0 (VIOLET) 35mm 1/2 CIRCLE TAPER POINT NEEDLE	DC0M35-120	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8681677075861	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 90cm USP SIZE 0 (VIOLET) 35mm 3/8 CIRCLE REVERSE CUTTING NEEDLE	DC0K35-90	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8681677075878	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 120cm USP SIZE 3-0 (VIOLET) 20mm 1/2 CIRCLE TAPER POINT NEEDLE	DC3M20-120	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8681677075977	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 100cm USP SIZE 2-0 (VIOLET) 35mm 3/8 CIRCLE REVERSE CUTTING NEEDLE	DC2K35-100	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8681677075984	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 75cm USP SIZE 5-0 (VIOLET) 18mm 3/8 CIRCLE REVERSE CUTTING NEEDLE	DC5K18	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8681677075991	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 75cm USP SIZE 2-0 (VIOLET) 20mm 3/8 CIRCLE REVERSE CUTTING NEEDLE	DC2K20	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60

	GMD Group Medikal San. Ve Tic. A.Ş		Ref No.	DC-06
	Department	Device master file	Rev No.	03
	Damacryl and PGA	PGA (Absorbable)	Issue Date	21.07.2023


8681677076165	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 75cm USP SIZE 6-0 (VIOLET) 16mm 1/2 CIRCLE TAPER POINT NEEDLE	DC6M16	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8681677076226	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 90cm USP SIZE 0 (VIOLET) 37mm 1/2 CIRCLE TAPER POINT NEEDLE	DC0M37-90	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8681677076257	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 3PCS X 45cm USP SIZE 1 (VIOLET) WITHOUT NEEDLE	DC1-3X45	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8681677076264	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 3PCS X 45cm USP SIZE 2 (VIOLET) WITHOUT NEEDLE	DC02-3X45	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8681677076271	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 6PCS X 45cm USP SIZE 1 (VIOLET) WITHOUT NEEDLE	DC1-6X45	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8681677076288	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 6PCS X 45cm USP SIZE 2 (VIOLET) WITHOUT NEEDLE	DC02-6X45	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8681677076639	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 5PCS X 45cm USP SIZE 0 (VIOLET) WITHOUT NEEDLE	DC0-5X45	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60

	GMD Group Medikal San. Ve Tic. A.Ş		Ref No.	DC-06
	Department	Device master file	Rev No.	03
	Damacryl and PGA	PGA (Absorbable)	Issue Date	21.07.2023


8681677076677	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 50cm USP SIZE 0 (VIOLET) WITHOUT NEEDLE	DC0-1X50	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8681677076646	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 5PCS X 75cm USP SIZE 0 (VIOLET) WITHOUT NEEDLE	DC0-5X75	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8681677076653	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 100cm USP SIZE 0 (VIOLET) 40mm 1/2 CIRCLE NORMAL CUTTING EDGE NEEDLE	DC0KDN40-100	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8681677076660	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 5PCS X 75cm USP SIZE 1 (VIOLET) WITHOUT NEEDLE	DC1-5X75	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8681677076684	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 75cm USP SIZE 1 (VIOLET) 36mm 1/2 CIRCLE TAPER POINT NEEDLE	DC1M36	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8681677076691	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 5PCS X 45cm USP SIZE 4-0 (VIOLET) WITHOUT NEEDLE	DC4-5X45	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8681677076707	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 90cm USP SIZE 4-0 (VIOLET) 22mm 1/2 CIRCLE TAPER POINT NEEDLE	DC4M22-90	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60

	GMD Group Medikal San. Ve Tic. A.Ş		Ref No.	DC-06
	Department	Device master file	Rev No.	03
	Damacryl and PGA	PGA (Absorbable)	Issue Date	21.07.2023


8681677076844	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 100cm USP SIZE 2-0 (VIOLET) 40mm 3/8 CIRCLE REVERSE CUTTING NEEDLE	DC2K40-100	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8681677076875	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 75cm USP SIZE 3-0 (VIOLET) 24mm 3/8 CIRCLE REVERSE CUTTING NEEDLE	DC3K24	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8681677076929	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 120cm USP SIZE 0 (VIOLET) 40mm 1/2 CIRCLE TAPER POINT NEEDLE	DC0M40-120	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8681677076998	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 75cm USP SIZE 3-0 (VIOLET) 17mm 1/2 CIRCLE TAPER POINT NEEDLE	DC3M17	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8681677077018	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 150cm USP SIZE 3-0 (VIOLET DYED) WITHOUT NEEDLE	DC3-15	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8681677077056	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 75cm USP SIZE 6-0 (VIOLET) 10mm 1/2 CIRCLE TAPER POINT NEEDLE	DC6M10	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8681677077162	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURES 1PCS X 180cm USP SIZE 1 (VIOLET) 50mm 1/2 CIRCLE TAPER POINT & 50mm 3/8 CIRCLE REVERSE CUTTING DOUBLE NEEDLE	DC1K50M50-D-180	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60

	GMD Group Medikal San. Ve Tic. A.Ş		Ref No.	DC-06
	Department	Device master file	Rev No.	03
	Damacryl and PGA	PGA (Absorbable)	Issue Date	21.07.2023


8681677077278	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 5PCS X 45cm USP SIZE 1 (VIOLET) WITHOUT NEEDLE	DC1-5X45	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8681677077285	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 75cm USP SIZE 2-0 (VIOLET) 24mm 3/8 CIRCLE REVERSE CUTTING NEEDLE	DC2K24	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8681677077292	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 100cm USP SIZE 2-0 (VIOLET) 35mm 1/2 CIRCLE TAPER POINT NEEDLE	DC2M35-100	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8681677077308	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 100cm USP SIZE 0 (VIOLET) WITHOUT NEEDLE	DC0-1X100	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8681677077353	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 75cm USP SIZE 3-0 (VIOLET) 35mm 3/8 CIRCLE REVERSE CUTTING NEEDLE	DC3K35	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8681677077391	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 75cm USP SIZE 0 (VIOLET) 37mm 1/2 CIRCLE TAPER POINT NEEDLE	DC0M37	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8681677077452	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 75cm USP SIZE 3-0 (VIOLET) 37mm 1/2 CIRCLE TAPER POINT NEEDLE	DC3M37	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60

	GMD Group Medikal San. Ve Tic. A.Ş		Ref No.	DC-06
	Department	Device master file	Rev No.	03
	Damacryl and PGA	PGA (Absorbable)	Issue Date	21.07.2023


8681677077742	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 75cm USP SIZE 2-0 (VIOLET) 45mm 1/2 CIRCLE TAPER POINT NEEDLE	DC2M45	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8681677077766	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 75cm USP SIZE 4-0 (VIOLET) 22mm 3/8 CIRCLE REVERSE CUTTING NEEDLE	DC4K22	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8681677077940	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 5PCS X 45cm USP SIZE 2 (VIOLET) WITHOUT NEEDLE	DC02-5X45	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8681677077957	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 75cm USP SIZE 5-0 (VIOLET) 19mm 3/8 CIRCLE REVERSE CUTTING NEEDLE	DC5K19	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8681677077964	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 75cm USP SIZE 5-0 (VIOLET) 17mm 1/2 CIRCLE TAPER POINT NEEDLE	DC5M17	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8681677078008	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 75cm USP SIZE 2-0 (VIOLET) 24mm 1/2 CIRCLE TAPER POINT NEEDLE	DC2M24	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8681677078176	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 100cm USP SIZE 3-0 (VIOLET) 30mm 3/8 CIRCLE REVERSE CUTTING NEEDLE	DC3K30-100	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60

	GMD Group Medikal San. Ve Tic. A.Ş		Ref No.	DC-06
	Department	Device master file	Rev No.	03
	Damacryl and PGA	PGA (Absorbable)	Issue Date	21.07.2023


8681677078190	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 75cm USP SIZE 3-0 (VIOLET) 20mm 3/8 CIRCLE REVERSE CUTTING NEEDLE	DC3K20	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8681677078183	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 75cm USP SIZE 4-0 (VIOLET) 15mm 1/2 CIRCLE TAPER POINT NEEDLE	DC4M15	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8681677078282	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 75cm USP SIZE 4-0 (VIOLET) 18mm 3/8 CIRCLE REVERSE CUTTING NEEDLE	DC4K18	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8681677078664	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 75cm USP SIZE 0 (VIOLET) 26mm 1/2 CIRCLE TAPER POINT NEEDLE	DC0HR26	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8681677078671	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 75cm USP SIZE 0 (VIOLET) 35mm 3/8 CIRCLE REVERSE CUTTING NEEDLE	DC0DS35	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8681677078688	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 90cm USP SIZE 2 (VIOLET) 45mm 1/2 CIRCLE TAPER POINT NEEDLE	DC02HR45-90	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8681677078909	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 75cm USP SIZE 6-0 (VIOLET) 16mm 3/8 CIRCLE REVERSE CUTTING NEEDLE	DC6DS16	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60

	GMD Group Medikal San. Ve Tic. A.Ş		Ref No.	DC-06
	Department	Device master file	Rev No.	03
	Damacryl and PGA	PGA (Absorbable)	Issue Date	21.07.2023


8681677078916	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 100cm USP SIZE 3-0 (VIOLET) 35mm 3/8 CIRCLE REVERSE CUTTING NEEDLE	DC3DS35-100	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8681677078923	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 100cm USP SIZE 3-0 (VIOLET) 35mm 1/2 CIRCLE TAPER POINT NEEDLE	DC3HR35-100	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8681677078930	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 75cm USP SIZE 5-0 (VIOLET) 20mm 3/8 CIRCLE REVERSE CUTTING NEEDLE	DC5DS20	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8681677078947	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 75cm USP SIZE 5-0 (VIOLET) 20mm 1/2 CIRCLE TAPER POINT NEEDLE	DC5HR20	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8681677079036	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 90cm USP SIZE 2 (VIOLET) 65mm 1/2 CIRCLE TAPER POINT NEEDLE	DC02HR65-90	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8681677079043	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 75cm USP SIZE 1 (VIOLET) 31mm J SHAPE TAPER CUT NEEDLE	DC1JSY31	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8681677079241	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 75cm USP SIZE 4-0 (VIOLET) 30mm 1/2 CIRCLE TAPER POINT NEEDLE	DC40HR30	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60

	GMD Group Medikal San. Ve Tic. A.Ş		Ref No.	DC-06
	Department	Device master file	Rev No.	03
	Damacryl and PGA	PGA (Absorbable)	Issue Date	21.07.2023


8681677079258	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 75cm USP SIZE 2 (VIOLET) 35mm 1/2 CIRCLE TAPER POINT NEEDLE	DC02HR35	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8681677079494	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 75cm USP SIZE 3-0 (VIOLET) 23mm 1/2 CIRCLE TAPER POINT NEEDLE	DC30HR23	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8681677079500	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 75cm USP SIZE 3-0 (VIOLET) 35mm 1/2 CIRCLE TAPER POINT NEEDLE	DC30HR35	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8681677079517	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 75cm USP SIZE 4-0 (VIOLET) 18mm 1/2 CIRCLE TAPER POINT NEEDLE	DC40HR18	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8681677079630	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 75cm USP SIZE 2-0 (UNDYED) 30mm 1/2 CIRCLE TAPER POINT NEEDLE	DC20HR30U	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8681677079647	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 150cm USP SIZE 2-0 (UNDYED) WITHOUT NEEDLE	DC20-150U	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8684154200022	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 75cm USP SIZE 2-0 (UNDYED) 24mm 1/2 CIRCLE TAPER POINT NEEDLE	DC20HR24U	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60

	GMD Group Medikal San. Ve Tic. A.Ş		Ref No.	DC-06
	Department	Device master file	Rev No.	03
	Damacryl and PGA	PGA (Absorbable)	Issue Date	21.07.2023


8681677079661	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 75cm USP SIZE 2-0 (UNDYED) 27mm 1/2 CIRCLE TAPER POINT NEEDLE	DC20HR27U	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8681677079678	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 150cm USP SIZE 3-0 (UNDYED) WITHOUT NEEDLE	DC30-150U	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8681677079685	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 75cm USP SIZE 3-0 (UNDYED) 26mm 1/2 CIRCLE TAPER POINT NEEDLE	DC30HR26U	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8681677079708	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 75cm USP SIZE 3-0 (VIOLET) 16mm 3/8 CIRCLE REVERSE CUTTING NEEDLE	DC30DS16	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8681677079715	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 75cm USP SIZE 4-0 (VIOLET) 17mm 3/8 CIRCLE REVERSE CUTTING NEEDLE	DC40DS17	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8684154200145	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 75cm USP SIZE 2-0 (VIOLET) 22mm 1/2 CIRCLE TAPER POINT NEEDLE	DC20HR22	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8684154200169	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 75cm USP SIZE 0 (VIOLET) 40mm 1/2 CIRCLE NORMAL CUTTING EDGE NEEDLE	DC00HC40	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60

	GMD Group Medikal San. Ve Tic. A.Ş		Ref No.	DC-06
	Department	Device master file	Rev No.	03
	Damacryl and PGA	PGA (Absorbable)	Issue Date	21.07.2023


8684154200176	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 75cm USP SIZE 3-0 (VIOLET) 25mm 3/8 CIRCLE REVERSE CUTTING NEEDLE	DC30DS25	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8684154200183	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 75cm USP SIZE 5-0 (VIOLET) 17mm 1/2 CIRCLE REVERSE CUTTING NEEDLE	DC50HS17	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8684154200213	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 75cm USP SIZE 6-0 (VIOLET) 16mm 1/2 CIRCLE REVERSE CUTTING NEEDLE	DC60HS16	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8684154200664	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 75cm USP SIZE 4-0 (VIOLET DYED) 16mm 3/8 CIRCLE REVERSE CUTTING NEEDLE	DC40DS16	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8684154200671	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 75cm USP SIZE 4-0 (VIOLET DYED) 19mm 3/8 CIRCLE REVERSE CUTTING NEEDLE	DC40DS19	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8684154200688	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 75cm USP SIZE 4-0 (VIOLET) 17mm 1/2 CIRCLE REVERSE CUTTING NEEDLE	DC40HS17	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8684154200695	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 75cm USP SIZE 4-0 (VIOLET) 20mm 1/2 CIRCLE REVERSE CUTTING NEEDLE	DC40HS20	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60

	GMD Group Medikal San. Ve Tic. A.Ş		Ref No.	DC-06
	Department	Device master file	Rev No.	03
	Damacryl and PGA	PGA (Absorbable)	Issue Date	21.07.2023


8684154200701	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 75cm USP SIZE 4-0 (VIOLET) 18mm 1/2 CIRCLE TAPER CUT NEEDLE	DC40HRT18	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8684154200718	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 75cm USP SIZE 4-0 (VIOLET DYED) 17mm 1/2 CIRCLE TAPER POINT NEEDLE	DC40HR17	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8684154200725	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 75cm USP SIZE 4-0 (VIOLET DYED) 20mm 1/2 CIRCLE TAPER POINT NEEDLE	DC40HR20	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8684154200732	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 75cm USP SIZE 4-0 (VIOLET) 25mm COMPOUND CURVE TAPER CUT NEEDLE	DC40ART25	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8684154200749	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 75cm USP SIZE 3-0 (VIOLET DYED) 19mm 3/8 CIRCLE REVERSE CUTTING NEEDLE	DC30DS19	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8684154200978	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 75cm USP SIZE 3-0 (VIOLET) 16mm 1/2 CIRCLE TAPER CUT NEEDLE	DC30HRT16	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8684154201036	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 45cm USP SIZE 4-0 (VIOLET) 16mm 3/8 CIRCLE NORMAL CUTTING NEEDLE	DC40DC16-45	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60

	GMD Group Medikal San. Ve Tic. A.Ş		Ref No.	DC-06
	Department	Device master file	Rev No.	03
	Damacryl and PGA	PGA (Absorbable)	Issue Date	21.07.2023


8684154201111	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 75cm USP SIZE 6-0 (VIOLET) 19mm 3/8 CIRCLE REVERSE CUTTING NEEDLE	DC60DS19	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8684154201128	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 7PCS X 75cm USP SIZE 4-0 (VIOLET) WITHOUT NEEDLE	DC40-7X75	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8684154201135	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 7PCS X 75cm USP SIZE 3-0 (VIOLET) WITHOUT NEEDLE	DC30-7X75	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8684154201500	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 75cm USP SIZE 1 (VIOLET) 26mm 1/2 CIRCLE TAPER POINT NEEDLE	DC01HR26	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8684154201517	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 75cm USP SIZE 2-0 (VIOLET) 30mm 1/2 CIRCLE REVERSE CUTTING NEEDLE	DC20HS30	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8684154201623	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 250cm USP SIZE 0 (VIOLET) WITHOUT NEEDLE	DC00-250	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8684154201630	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 250cm USP SIZE 1 (VIOLET) WITHOUT NEEDLE	DC01-250	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60

	GMD Group Medikal San. Ve Tic. A.Ş		Ref No.	DC-06
	Department	Device master file	Rev No.	03
	Damacryl and PGA	PGA (Absorbable)	Issue Date	21.07.2023


8684154201647	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 250cm USP SIZE 2 (VIOLET) WITHOUT NEEDLE	DC02-250	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8684154201692	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 75cm USP SIZE 3-0 (VIOLET) 26mm 1/2 CIRCLE REVERSE CUTTING NEEDLE	DC30HS26	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8684154201708	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 75cm USP SIZE 4-0 (VIOLET) 26mm 1/2 CIRCLE REVERSE CUTTING NEEDLE	DC40HS26	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8684154201760	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 45cm USP SIZE 3-0 (VIOLET) 16mm 3/8 CIRCLE NORMAL CUTTING NEEDLE	DC30DC16-45	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8684154202811	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 75cm USP SIZE 2-0 (VIOLET) 19mm 3/8 CIRCLE REVERSE CUTTING NEEDLE	DC20DS19	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8684154202828	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 75cm USP SIZE 2-0 (VIOLET) WITHOUT NEEDLE	DC20-1X75	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8684154202880	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 90cm USP SIZE 0 (VIOLET) 40mm 1/2 CIRCLE REVERSE CUTTING NEEDLE	DC00HS40-90	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60

	GMD Group Medikal San. Ve Tic. A.Ş		Ref No.	DC-06
	Department	Device master file	Rev No.	03
	Damacryl and PGA	PGA (Absorbable)	Issue Date	21.07.2023


8684154202897	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 100cm USP SIZE 1 (VIOLET) 40mm 1/2 CIRCLE REVERSE CUTTING NEEDLE	DC01HS40-100	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8684154202903	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 100cm USP SIZE 1 (VIOLET) 40mm 3/8 CIRCLE REVERSE CUTTING NEEDLE	DC01DS40-100	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8684154202910	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 90cm USP SIZE 2-0 (VIOLET) 30mm 1/2 CIRCLE REVERSE CUTTING NEEDLE	DC20HS30-90	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8684154202927	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 100cm USP SIZE 2 (VIOLET) 40mm 3/8 CIRCLE REVERSE CUTTING NEEDLE	DC02DS40-100	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8684154203269	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 75cm USP SIZE 6-0 (VIOLET) 19mm 1/2 CIRCLE TAPER POINT NEEDLE	DC60HR19	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8684154203276	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 75cm USP SIZE 2-0 (VIOLET) 36mm 3/8 CIRCLE REVERSE CUTTING NEEDLE	DC20DS36	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8684154203436	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 45cm USP SIZE 6-0 (VIOLET) 13mm 3/8 CIRCLE REVERSE CUTTING NEEDLE	DC60DS13-45	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60

	GMD Group Medikal San. Ve Tic. A.Ş		Ref No.	DC-06
	Department	Device master file	Rev No.	03
	Damacryl and PGA	PGA (Absorbable)	Issue Date	21.07.2023


8684154203825	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 75cm USP SIZE 3-0 (VIOLET) 18mm 3/8 CIRCLE REVERSE CUTTING NEEDLE	DC30DS18	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8684154203832	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 75cm USP SIZE 4-0 (VIOLET) 24mm 3/8 CIRCLE REVERSE CUTTING NEEDLE	DC40DS24	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8684154203849	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 75cm USP SIZE 5-0 (VIOLET) 24mm 3/8 CIRCLE REVERSE CUTTING NEEDLE	DC50DS24	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8684154203856	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 120cm USP SIZE 0 (VIOLET) 40mm 1/2 CIRCLE REVERSE CUTTING NEEDLE	DC00HS40-120	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8684154203863	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 120cm USP SIZE 2 (VIOLET) 40mm 1/2 CIRCLE REVERSE CUTTING NEEDLE	DC02HS40-120	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8684154203870	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 120cm USP SIZE 2 (VIOLET) 48mm 1/2 CIRCLE REVERSE CUTTING NEEDLE	DC02HS48-120	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8684154203887	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 150cm USP SIZE 2 (VIOLET) 40mm 1/2 CIRCLE REVERSE CUTTING NEEDLE	DC02HS40-150	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60

	GMD Group Medikal San. Ve Tic. A.Ş		Ref No.	DC-06
	Department	Device master file	Rev No.	03
	Damacryl and PGA	PGA (Absorbable)	Issue Date	21.07.2023


8684154203894	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 120cm USP SIZE 1 (VIOLET) 30mm 3/8 CIRCLE REVERSE CUTTING NEEDLE	DC01DS30-120	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8684154203993	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 6PCS X 45cm USP SIZE 3-0 (VIOLET) WITHOUT NEEDLE	DC30-6X45	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8684154204006	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 45cm USP SIZE 3-0 (VIOLET) 18mm 1/2 CIRCLE TAPER POINT NEEDLE	DC30HR18-45	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8684154204013	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 45cm USP SIZE 4-0 (VIOLET) 18mm 1/2 CIRCLE TAPER POINT NEEDLE	DC40HR18-45	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8684154204020	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 45cm USP SIZE 5-0 (VIOLET) 19mm 3/8 CIRCLE REVERSE CUTTING NEEDLE	DC50DS19-45	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8684154204259	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 75cm USP SIZE 2 (VIOLET) 65mm 1/2 CIRCLE TAPER POINT NEEDLE	DC02HR65	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8681677085815	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 150cm USP SIZE 3-0 (VIOLET) 20mm 1/2 CIRCLE TAPER POINT NEEDLE	DC30HR20-150	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60

	GMD Group Medikal San. Ve Tic. A.Ş		Ref No.	DC-06
	Department	Device master file	Rev No.	03
	Damacryl and PGA	PGA (Absorbable)	Issue Date	21.07.2023


8684154204358	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 120cm USP SIZE 2 (VIOLET) 30mm 1/2 CIRCLE TAPER POINT NEEDLE	DC02HR30-120	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8684154204365	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 120cm USP SIZE 2 (VIOLET) 30mm 3/8 CIRCLE REVERSE CUTTING NEEDLE	DC02DS30-120	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8684154204372	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 120cm USP SIZE 2 (VIOLET) 50mm 1/2 CIRCLE REVERSE CUTTING NEEDLE	DC02HS50-120	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8684154204389	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 120cm USP SIZE 1 (VIOLET) 40mm 1/2 CIRCLE REVERSE CUTTING NEEDLE	DC01HS40-120	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8684154204396	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 75cm USP SIZE 2-0 (VIOLET) 30mm 3/8 CIRCLE TAPER POINT NEEDLE	DC20DR30	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8684154204402	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 75cm USP SIZE 3-0 (VIOLET) 20mm 3/8 CIRCLE TAPER POINT NEEDLE	DC30DR20	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8684154204495	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 150cm USP SIZE 3 (VIOLET) WITHOUT NEEDLE	DC03-1X150	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60

	GMD Group Medikal San. Ve Tic. A.Ş		Ref No.	DC-06
	Department	Device master file	Rev No.	03
	Damacryl and PGA	PGA (Absorbable)	Issue Date	21.07.2023


8684154204624	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 75cm USP SIZE 4-0 (VIOLET) 19mm 1/2 CIRCLE REVERSE CUTTING NEEDLE	DC40HS19	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8684154204648	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 150cm USP SIZE 2 (VIOLET) 40mm 3/8 CIRCLE REVERSE CUTTING NEEDLE	DC02DS40-150	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8684154204709	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 6PCS X 45cm USP SIZE 2-0 (VIOLET) WITHOUT NEEDLE	DC20-6X45	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8684154204716	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 6PCS X 45cm USP SIZE 4-0 (VIOLET) WITHOUT NEEDLE	DC40-6X45	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8684154204723	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 5PCS X 50cm USP SIZE 2-0 (VIOLET) WITHOUT NEEDLE	DC20-5X50	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8684154206017	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 75cm USP SIZE 0 (VIOLET) 61mm 3/8 CIRCLE ROUND BODY (BLUNT TIP) NEEDLE	DC00DRN61	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8684154206024	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 75cm USP SIZE 1 (VIOLET) 61mm 3/8 CIRCLE ROUND BODY (BLUNT TIP) NEEDLE	DC01DRN61	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60

	GMD Group Medikal San. Ve Tic. A.Ş		Ref No.	DC-06
	Department	Device master file	Rev No.	03
	Damacryl and PGA	PGA (Absorbable)	Issue Date	21.07.2023


8684154206031	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 75cm USP SIZE 2 (VIOLET) 65mm 3/8 CIRCLE ROUND BODY (BLUNT TIP) NEEDLE	DC02DRN65	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8684154206055	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 10PCS X 75cm USP SIZE 2-0 (VIOLET) WITHOUT NEEDLE	DC20-10X75	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8684154206062	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 10PCS X 75cm USP SIZE 3-0 (VIOLET) WITHOUT NEEDLE	DC30-10X75	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8684154206079	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 75cm USP SIZE 3-0 (VIOLET) 24mm 1/2 CIRCLE TAPER POINT NEEDLE	DC30HR24	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8684154206086	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 75cm USP SIZE 4-0 (VIOLET) 19mm 1/2 CIRCLE TAPER POINT NEEDLE	DC40HR19	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8684154206185	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 45cm USP SIZE 6-0 (VIOLET) 8mm 1/4 CIRCLE SPATULA MICRO POINT DOUBLE NEEDLE	DC60VSPM8-D-45	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8684154206192	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 45cm USP SIZE 7-0 (VIOLET) 7mm 3/8 CIRCLE SPATULA MICRO POINT DOUBLE NEEDLE	DC70DSPM7-D-45	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60

	GMD Group Medikal San. Ve Tic. A.Ş		Ref No.	DC-06
	Department	Device master file	Rev No.	03
	Damacryl and PGA	PGA (Absorbable)	Issue Date	21.07.2023


8684154206437	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 90cm USP SIZE 2 (VIOLET) 40mm 1/2 CIRCLE REVERSE CUTTING NEEDLE	DC02HS40-90	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8684154206444	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 90cm USP SIZE 1 (VIOLET) 40mm 1/2 CIRCLE REVERSE CUTTING NEEDLE	DC01HS40-90	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8684154206840	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 90cm USP SIZE 2-0 (VIOLET) 26mm 3/8 CIRCLE REVERSE CUTTING NEEDLE	DC20DS26-90	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8684154207250	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 90cm USP SIZE 4-0 (VIOLET) 24mm 3/8 CIRCLE REVERSE CUTTING NEEDLE	DC40DS24-90	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8684154207267	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 90cm USP SIZE 3-0 (VIOLET) 24mm 3/8 CIRCLE REVERSE CUTTING NEEDLE	DC30DS24-90	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8684154207274	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 90cm USP SIZE 2-0 (VIOLET) 24mm 3/8 CIRCLE REVERSE CUTTING NEEDLE	DC20DS24-90	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8684154207281	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 90cm USP SIZE 2 (VIOLET) 30mm 3/8 CIRCLE REVERSE CUTTING NEEDLE	DC02DS30-90	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60

	GMD Group Medikal San. Ve Tic. A.Ş		Ref No.	DC-06
	Department	Device master file	Rev No.	03
	Damacryl and PGA	PGA (Absorbable)	Issue Date	21.07.2023


8684154207298	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 90cm USP SIZE 2 (VIOLET) 50mm 3/8 CIRCLE REVERSE CUTTING NEEDLE	DC02DS50-90	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8684154207304	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 90cm USP SIZE 1 (VIOLET) 31mm J SHAPE TAPER CUT NEEDLE	DC01JSY31-90	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8684154207311	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 75cm USP SIZE 2-0 (VIOLET) 25mm 1/2 CIRCLE TAPER POINT NEEDLE	DC20HR25	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8684154207328	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 300cm USP SIZE 1 (VIOLET) WITHOUT NEEDLE	DC01-1X300	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8684154207335	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 300cm USP SIZE 2 (VIOLET) WITHOUT NEEDLE	DC02-1X300	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8684154207342	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 100cm USP SIZE 0 (VIOLET) 30mm 1/2 CIRCLE TAPER POINT NEEDLE	DC00HR30-100	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8684154207403	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 75cm USP SIZE 3 (VIOLET) 40mm 1/2 CIRCLE TAPER POINT NEEDLE	DC03HR40	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60

	GMD Group Medikal San. Ve Tic. A.Ş		Ref No.	DC-06
	Department	Device master file	Rev No.	03
	Damacryl and PGA	PGA (Absorbable)	Issue Date	21.07.2023


8684154207410	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 75cm USP SIZE 3 (VIOLET) 30mm 1/2 CIRCLE TAPER POINT NEEDLE	DC03HR30	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8684154207427	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 75cm USP SIZE 2 (VIOLET) 45mm 1/2 CIRCLE TAPER POINT (HEAVY) NEEDLE	DC02HRH45	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8684154207434	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 75cm USP SIZE 1 (VIOLET) 37mm 1/2 CIRCLE TAPER POINT NEEDLE	DC01HR37	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8684154207441	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 75cm USP SIZE 0 (VIOLET) 40mm 1/2 CIRCLE REVERSE CUTTING NEEDLE	DC00HS40	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8684154207458	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 75cm USP SIZE 2-0 (VIOLET) 26mm 1/2 CIRCLE REVERSE CUTTING NEEDLE	DC20HS26	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8684154207465	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 75cm USP SIZE 4-0 (VIOLET) 37mm 1/2 CIRCLE TAPER POINT NEEDLE	DC40HR37	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8684154207908	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 45cm USP SIZE 6-0 (VIOLET) 10mm 3/8 CIRCLE REVERSE CUTTING NEEDLE	DC60DS10-45	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60

	GMD Group Medikal San. Ve Tic. A.Ş		Ref No.	DC-06
	Department	Device master file	Rev No.	03
	Damacryl and PGA	PGA (Absorbable)	Issue Date	21.07.2023


8684154208097	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 90cm USP SIZE 4-0 (VIOLET) 26mm 1/2 CIRCLE TAPER POINT NEEDLE	DC40HR26-90	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8684154208103	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 90cm USP SIZE 5-0 (VIOLET) 26mm 1/2 CIRCLE TAPER POINT NEEDLE	DC50HR26-90	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8684154208233	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 75cm USP SIZE 0 (VIOLET) 26mm 3/8 CIRCLE REVERSE CUTTING NEEDLE	DC00DS26	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8684154208240	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 75cm USP SIZE 0 (VIOLET) 26mm 1/2 CIRCLE REVERSE CUTTING NEEDLE	DC00HS26	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8684154208257	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 75cm USP SIZE 4-0 (VIOLET) 16mm 1/2 CIRCLE REVERSE CUTTING NEEDLE	DC40HS16	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8684154208264	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 75cm USP SIZE 3-0 (VIOLET) 16mm 1/2 CIRCLE TAPER POINT NEEDLE	DC30HR16	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8684154208271	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 75cm USP SIZE 3-0 (VIOLET) 20mm 1/2 CIRCLE REVERSE CUTTING NEEDLE	DC30HS20	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60

	GMD Group Medikal San. Ve Tic. A.Ş		Ref No.	DC-06
	Department	Device master file	Rev No.	03
	Damacryl and PGA	PGA (Absorbable)	Issue Date	21.07.2023

8684154208288	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 75cm USP SIZE 3-0 (VIOLET) 16mm 1/2 CIRCLE REVERSE CUTTING NEEDLE	DC30HS16	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8684154208295	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 75cm USP SIZE 5-0 (VIOLET) 16mm 1/2 CIRCLE REVERSE CUTTING NEEDLE	DC50HS16	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8684154208301	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 75cm USP SIZE 5-0 (VIOLET) 20mm 1/2 CIRCLE REVERSE CUTTING NEEDLE	DC50HS20	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8684154208400	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 75cm USP SIZE 3 (VIOLET) 48mm 1/2 CIRCLE TAPER POINT NEEDLE	DC03HR48	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8684154208554	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 75cm USP SIZE 3-0 (VIOLET) 26mm 1/2 CIRCLE TAPER POINT DOUBLE NEEDLE	DC30HR26-D	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8684154208868	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 75cm USP SIZE 5-0 (VIOLET) 19mm 3/8 CIRCLE REVERSE CUTTING (PRIME POINT) NEEDLE	DC50DSPR19	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8684154208875	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 75cm USP SIZE 4-0 (VIOLET) 19mm 3/8 CIRCLE REVERSE CUTTING (PRIME POINT) NEEDLE	DC40DSPR19	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60

	GMD Group Medikal San. Ve Tic. A.Ş		Ref No.	DC-06
	Department	Device master file	Rev No.	03
	Damacryl and PGA	PGA (Absorbable)	Issue Date	21.07.2023

8684154208882	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 75cm USP SIZE 3-0 (VIOLET) 25mm 3/8 CIRCLE REVERSE CUTTING (PRIME POINT) NEEDLE	DC30DSPR25	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8684154208899	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 300cm USP SIZE 0 (VIOLET) WITHOUT NEEDLE	DC00-1X300	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8684154208950	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 250cm USP SIZE 3-0 (VIOLET) WITHOUT NEEDLE	DC30-1X250	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8684154214739	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 75cm USP SIZE 0 (VIOLET) 26mm 1/2 CIRCLE NORMAL CUTTING EDGE NEEDLE	DC00HC26	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8684154214760	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 75cm USP SIZE 2-0 (VIOLET) 26mm 1/2 CIRCLE NORMAL CUTTING EDGE NEEDLE	DC20HC26	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8684154214777	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 75cm USP SIZE 3-0 (VIOLET) 16mm 1/2 CIRCLE NORMAL CUTTING EDGE NEEDLE	DC30HC16	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
8684154214784	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 75cm USP SIZE 3-0 (VIOLET) 20mm 1/2 CIRCLE NORMAL CUTTING EDGE NEEDLE	DC30HC20	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60

	GMD Group Medikal San. Ve Tic. A.Ş		Ref No.	DC-06
	Department	Device master file	Rev No.	03
	Damacryl and PGA	PGA (Absorbable)	Issue Date	21.07.2023

8684154214791	DAMACRYL: PGA BRAIDED ABSORBABLE SUTURE 1PCS X 75cm USP SIZE 5-0 (VIOLET) 20mm 1/2 CIRCLE NORMAL CUTTING EDGE NEEDLE	DC50HC20	DAMACRYL	13908 - Suture, polyglycolic acid	Class III	60
---------------	--	----------	----------	-----------------------------------	-----------	----


Product specifications

POLY GLYCOLIC ACID (DAMACRYL)

EP SIZE	USP SIZE	DIA METER			Knot pull tensile strength		Suture-Needle tensile strength	
		A Dia meter average (mm)	B		C Average values (kgf)	D Min values (kgf)	A Average (kgf)	B Min values (kgf)
			Min values (mm)	Max values (mm)				
0.7	6-0	0.070-0.099	0.060	0.125	0.25	0.071	0.17	0.08
1	5-0	0.100-0.149	0.085	0.175	0.68	0.183	0.23	0.11
1.5	4-0	0.150-0.199	0.125	0.225	0.95	0.387	0.46	0.23
2	3-0	0.200-0.249	0.175	0.275	1.77	0.765	0.69	0.35
2.5	-	0.250-0.299	0.225	0.325	2.68	1.05	0.92	0.46
3	2-0	0.300-0.349	0.275	0.375	2.68	1.019	1.12	0.46
3.5	0	0.350-0.399	0.325	0.450	3.90	1.274	1.53	0.46
4	1	0.400-0.499	0.375	0.550	5.08	2.039	1.80	0.60
5	2	0.500-0.599	0.450	0.650	6.35	2.803	1.80	0.70

1.4.1 Components and Raw Materials of product: Surgical Suture Surgical Needle

Manufacturer/Supplier	Material
METABIOMED	Polyglycolic Acid (pga), beige (undyed), coated

	GMD Group Medikal San. Ve Tic. A.Ş		Ref No.	DC-06
	Department	Device master file	Rev No.	03
	Damacryl and PGA	PGA (Absorbable)	Issue Date	21.07.2023

Sima Ambalaj	Shrinking Cover
Quality Needles Pvt	Stainless Steel
Kurşun Ambalaj FTM Ambalaj	First packing (aluminum foil) Second packing (medical paper)
Endüstri Matbaacılık	Paper (slat)

1.5.Accessories Information+

There is no accessories.

1.7.Suitability of Raw Material

Raw materials are supplied from the suppliers given in Components of product. Raw materials analysis reports are given in Test Reports.

1.8.Indication, Contraindication and Adverse Reactions

INDICATIONS

DAMACRYL AND PGA suture is indicated for use in general soft tissue approximation and ligation, including the use in ophthalmic procedures. DAMACRYL AND PGA is not for use in cardiovascular or neurological tissues.

CONTRAINDICATIONS

DAMACRYL AND PGA suture, being absorbable, should not be used where extended approximation of tissue is required.

ADVERSE REACTIONS


Adverse effects associated with the use of this device include wound dehiscence, failure to provide adequate wound support in closure of the sites where expansion, stretching, or distension occur, failure to provide adequate wound support in elderly, malnourished or debilitated patients or in patients suffering from conditions which may delay wound healing, infection, minimal acute inflammatory tissue reaction, localized irritation when skin sutures are left in place for greater than 7 days, suture extrusion and delayed absorption in tissue with poor blood supply, calculi formation in urinary and biliary tracts when prolonged contact with salt solutions such as urine and bile occurs, and transitory local irritation at the wound site.

Broken needles may result in extended or additional surgeries or residual foreign bodies. Inadvertent needle sticks with contaminated surgical needles may result in the transmission of blood borne pathogens.

1.9. Material Specifications:

The following components and compounds are used for the manufacturing of the surgical sutures:

1.9.1.Thread (PGA multifilament, Mepfil-Lac ®)


	GMD Group Medikal San. Ve Tic. A.Ş		Ref No.	DC-06
	Department	Device master file	Rev No.	03
	Damacryl and PGA	PGA (Absorbable)	Issue Date	21.07.2023

The following specifications are met:

Appearance:

Beige(undyed) colored multifilament braided multifibers according to the size of the strands that meets all the requirements of the European Pharmacopoeia (EP) for Synthetic Absorbable Multifilament Sutures and the requirements of United States Pharmacopoeia (USP) Absorbable Sutures except for a slight oversize in diameter.

Formula of polyglycolic acid:

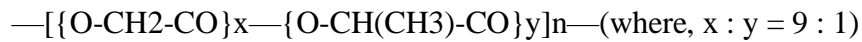

Content of residues:

Monomers ≤ 1.0 w/w % Catalysts \leq

100 ppm Heavy metals ≤ 10 ppm

Solvents ≤ 500 ppm

Additives:

Stabilizing agent: none Antioxidizing agent:

none **Dyestuffs**

Undyed braids: none

Coating

Polycaprolactone and calcium stearate




Content : ≤ 5 w/w %

Physical Specification

Damacryl and PGA braids meet all requirements established by the United States Pharmacopoeia (USP) and the European Pharmacopoeia (EP) except for diameter.

The diameter could be slightly over the specification.

	GMD Group Medikal San. Ve Tic. A.Ş		Ref No.	DC-06
	Department	Device master file	Rev No.	03
	Damacryl and PGA	PGA (Absorbable)	Issue Date	21.07.2023

Knot pull tensile strength linear pull tensile strength for non-sterile braid

USP 29 NF24 Page 2050

USP Size	Metric Size (gauge no.)	Limits on Average Diameter (mm)		Limits on Average Knot-Pull (except where otherwise specified) ^a Tensile Strength (in kgf) ^b			Limits on Average Knot-Pull (except where otherwise specified) ^a Tensile Strength (in N) ^b		
		Min.	Max.	Class I Min.	Class II Min.	Class III Min.	Class I Min.	Class II Min.	Class III Min.
12-0	0.01	0.001	0.009	0.001 ^a	—	0.002 ^a	0.01 ^a	—	0.02 ^a
11-0	0.1	0.010	0.019	0.006 ^a	0.005 ^a	0.02 ^a	0.06 ^a	0.05 ^a	0.20 ^a
10-0	0.2	0.020	0.029	0.019 ^a	0.014 ^a	0.06 ^a	0.194 ^a	0.14 ^a	0.59 ^a
9-0	0.3	0.030	0.039	0.043 ^a	0.029 ^a	0.07 ^a	0.424 ^a	0.28 ^a	0.68 ^a
8-0	0.4	0.040	0.049	0.06	0.04	0.11	0.59	0.39	1.08
7-0	0.5	0.050	0.069	0.11	0.06	0.16	1.08	0.59	1.57
6-0	0.7	0.070	0.099	0.20	0.11	0.27	1.96	1.08	2.65
5-0	1	0.10	0.149	0.40	0.23	0.54	3.92	2.26	5.30
4-0	1.5	0.15	0.199	0.60	0.46	0.82	5.88	4.51	8.04
3-0	2	0.20	0.249	0.96	0.66	1.36	9.41	6.47	13.3
2-0	3	0.30	0.339	1.44	1.02	1.80	14.1	10.0	17.6
0	3.5	0.35	0.399	2.16	1.45	3.40 ^a	21.2	14.2	33.3 ^a
1	4	0.40	0.499	2.72	1.81	4.76 ^a	26.7	17.8	46.7 ^a
2	5	0.50	0.599	3.52	2.54	5.90 ^a	34.5	24.9	57.8 ^a
3 and 4	6	0.60	0.699	4.88	3.68	9.11 ^a	47.8	36.1	89.3 ^a
5	7	0.70	0.799	6.16	—	11.4 ^a	60.4	—	112 ^a
6	8	0.80	0.899	7.28	—	13.6 ^a	71.4	—	133 ^a
7	9	0.90	0.999	9.04	—	15.9 ^a	88.6	—	156 ^a
8	10	1.00	1.099	—	—	18.2 ^a	—	—	178 ^a
9	11	1.100	1.199	—	—	20.5 ^a	—	—	201 ^a
10	12	1.200	1.299	—	—	22.8 ^a	—	—	224 ^a

^a The tensile strength of sizes smaller than USP size 8-0 (metric size 0.4) is measured by straight pull. The tensile strength of sizes larger than USP size 2-0 (metric size 3) of monofilament Class III (metallic) Nonabsorbable Surgical Suture is measured by straight pull.

Silver wire meets the tensile strength values of Class I Sutures but is tested in the same manner as Class III Sutures.


^b The limits on knot-pull tensile strength apply to Nonabsorbable Surgical Suture that has been sterilized. For nonsterile Sutures of Class I and Class II, the limits are 25% higher.

Table 1. Standard Needle Attachment for Absorbable and Nonabsorbable Sutures

Metric Size (Gauge No.)		Limits on Needle Attachment				
Absorbable Collagen Suture	Nonabsorbable and Synthetic Absorbable Sutures	USP Size	Average (in kgf) (Min.)	Individual (in kgf) (Min.)	Average (in N) (Min.)	Individual (in N) (Min.)
	0.1	11-0	0.007	0.005	0.069	0.049
	0.2	10-0	0.014	0.010	0.137	0.098
0.4	0.3	9-0	0.021	0.015	0.206	0.147
0.5	0.4	8-0	0.050	0.025	0.490	0.245
0.7	0.5	7-0	0.080	0.040	0.784	0.392
1	0.7	6-0	0.17	0.08	1.67	0.784
1.5	1	5-0	0.23	0.11	2.25	1.08
2	1.5	4-0	0.45	0.23	4.41	2.25
3	2	3-0	0.68	0.34	6.67	3.33
3.5	3	2-0	1.10	0.45	10.8	4.41
4	3.5	0	1.50	0.45	14.7	4.41
5	4	1	1.80	0.60	17.6	5.88
6 and larger	5 and larger	2 and larger	1.80	0.70	17.6	6.86

Moisture content: ≤ 400 ppm

Retention ratio of tensile strength

	GMD Group Medikal San. Ve Tic. A.Ş		Ref No.	DC-06
	Department	Device master file	Rev No.	03
	Damacryl and PGA	PGA (Absorbable)	Issue Date	21.07.2023

About 75% of the initial strength after two week.

1.9.2.NEEDLE:

Needle material is made from stainless steel and meets the internationally accepted specifications of surgically used metals.

1.9.3. Specification of primary and secondary packaging material

Primary packaging material is made of Aluminium foil.

Secondary packaging material is made from a pouch made of medical paper.

1.10. Patient Populations

The patients are used for all people who are need to ligate, close and cover suturable wounds under sterile conditions by initially suturing the wound due to gastrointestinal, obstetric and gynecologic and urologic surgery.

1.11. Intended User

The device should used by physician and/or surgeons.

1.12. Repeat Applications, Including Any Restrictions as to the Number or Duration of Reapplications


After surgery, the closing and covering can be repeated if an undesirable case would be occurred and/or the wound would open. There is not any restriction as to the number or duration of reapplications.

1.13. Single Use / Reusable

The device is only single-used.

Repeat Applications, Including Any Restrictions as to the Number or Duration of Reapplications

After surgery, the closing and covering can be repeated if an undesirable case would be occurred and/or the wound would open. There is not any restriction as to the number or duration of reapplications

	GMD Group Medikal San. Ve Tic. A.Ş		Ref No.	DC-06
	Department	Device master file	Rev No.	03
	Damacryl and PGA	PGA (Absorbable)	Issue Date	21.07.2023

1.14. Mechanism of Action / Principles of Operation

DAMACRYL AND PGA sutures elicit minimal acute inflammatory reaction in tissues, which is followed by gradual encapsulation of the suture by fibrous connective tissue. Progressive loss of tensile strength and eventual absorption of *DAMACRYL AND PGA* sterile, synthetic, absorbable suture occurs by means of hydrolysis, where the polymer degrades to glycolic acid which is subsequently absorbed and metabolized by the body. Absorption begins as a loss of tensile strength without appreciable loss of mass. The material holds approximately 65% of its tensile strength by the end of two weeks and approximately 45% by the third week. The absorption is completed in approximately from 60 to 90 days.

1.15. Biodegradation mechanism of the device

The literature describing the biodegradation mechanism for PGA polyglycolic acid is indicated in the annexes.(Annex-1 , Annex-2, Annex-3, Annex-4)

1.16. Biological and Toxicological Evaluation

The following section contains information on the preclinical evaluation of the raw materials and the toxicological and biocompatibility tests on the sterile finished product Damacryl.(TF13 Clinical and Biological Evaluation)


1.17. Biocompatibility Tests of the Raw Materials for Medical Device Damacryl and PGA(ISO10993-1:2009)

Specific safety evaluation programs follow International Organization for Standardization (ISO) 10993 standards and Food and Drug Administration (FDA) guidance (May 1, 1995). The table below is based on ISO 10993-1 Evaluation and testing, 2009 edition. The crosses (X) indicate data endpoints that can be necessary for a biological safety evaluation per ISO 10993-1, based on a risk analysis. Where existing data are adequate, additional testing is not required. (Clinical and Biological Evaluation)

1.18. Manufacturing Process:

Manufacturing is performed at the İkitelli OSB, İsdök Sosyal Tesis 2 Sk. B Blok 1 B 7 Başakşehir/İstanbul, Turkey. The manufacturing is made under clean room conditions. The productions consist of manually insertion of the thread into the needles. After testing products are wound in a paper, are packed into aluminium paper, then packed again in a medical paper and the products are ETO sterilized in-house.

The following major manufacturing steps are used:

	GMD Group Medikal San. Ve Tic. A.Ş		Ref No.	DC-06
	Department	Device master file	Rev No.	03
	Damacryl and PGA	PGA (Absorbable)	Issue Date	21.07.2023

Manufacturing Methods

- a) Threads are applied with tipping solution
- b) Threads are cut to the required length.
- c) Threads are manually attached to needles using purpose made attaching machines in ISO 8 Clean Room.

Packaging Procedures

- a) After attaching, the suture is reverse wound into a retaining pouch leaving the needle accessible for easy withdrawal of the thread.
- b) The wound sutures are packed in an aluminum foil inner envelope (primary packaging) and labeled in ISO 8 Clean Room, but not sealed at this stage.
- c) This inner aluminum envelope is then placed in an outer envelope (secondary packaging), which is sealed.

One package unit of Damacryl and PGA contains 12 single pouches.

Sterilization


Damacryl and PGA(PGA) sutures are sterilized in-house by ETO. They are then vacuum- dried, sealed and packed.

Sterilization process is verified by chemical indicators and biological indicators are also used routinely to challenge the procedure. The sterilization process is validated according to EN ISO 11135-1:2014. Main routine sterilization parameters are minimum 30 °C 30% Humidity with a gas injection period of 5 hours.

As regards the manufacturing / assembly area, environmental conditions are monitored at regular intervals using standard microbiological methods (i.e. culture plates).

Drying

Damacryl and PGA sutures are dried after ETO sterilization. They are vacuum-dried using a vacuum oven (3 hours, 60°C) Than with Nitrogen washing for 12 hours and a final vacuum 5 hours and drying for another 3 hours. At the end of the drying process the aluminium packages are sealed under controlled conditions. The sutures are then boxed, 12 per box and boxes are wrapped by shrinking cover.

	GMD Group Medikal San. Ve Tic. A.Ş		Ref No.	DC-06
	Department	Device master file	Rev No.	03
	Damacryl and PGA	PGA (Absorbable)	Issue Date	21.07.2023

Packaging into cartons

Twelve sutures pouches are packed into a carton and labelled. Each carton is polyolefin wrapped. Hundred cartons are packed in a shipper, which is sealed. Each shipper is allocated a unique sequential identification number and labelled with the contents and their batch numbers.

Production Hygiene:


Production (Assembly) is performed under ISO 8 clean room conditions. Bioburden is checked on a routinely basis and controlled by responsible persons.

ETO-Residual Tests

ISO 10993-7:2008 provides limits for sterilization residuals based on product application.

ETO residuals are tested according to the requirements of ISO 10993-7:2008 by using gas chromatography.

The residual test is carried out according to EN 10993-7 and the results are shown in Test Reports.


	GMD Group Medikal San. Ve Tic. A.Ş		Ref No.	DC-06
	Department	Device master file	Rev No.	03
	Damacryl and PGA	PGA (Absorbable)	Issue Date	21.07.2023


Manufacturing Flow Chart


WORKFLOW


RELATED DOCUMENT

SUPERVISOR

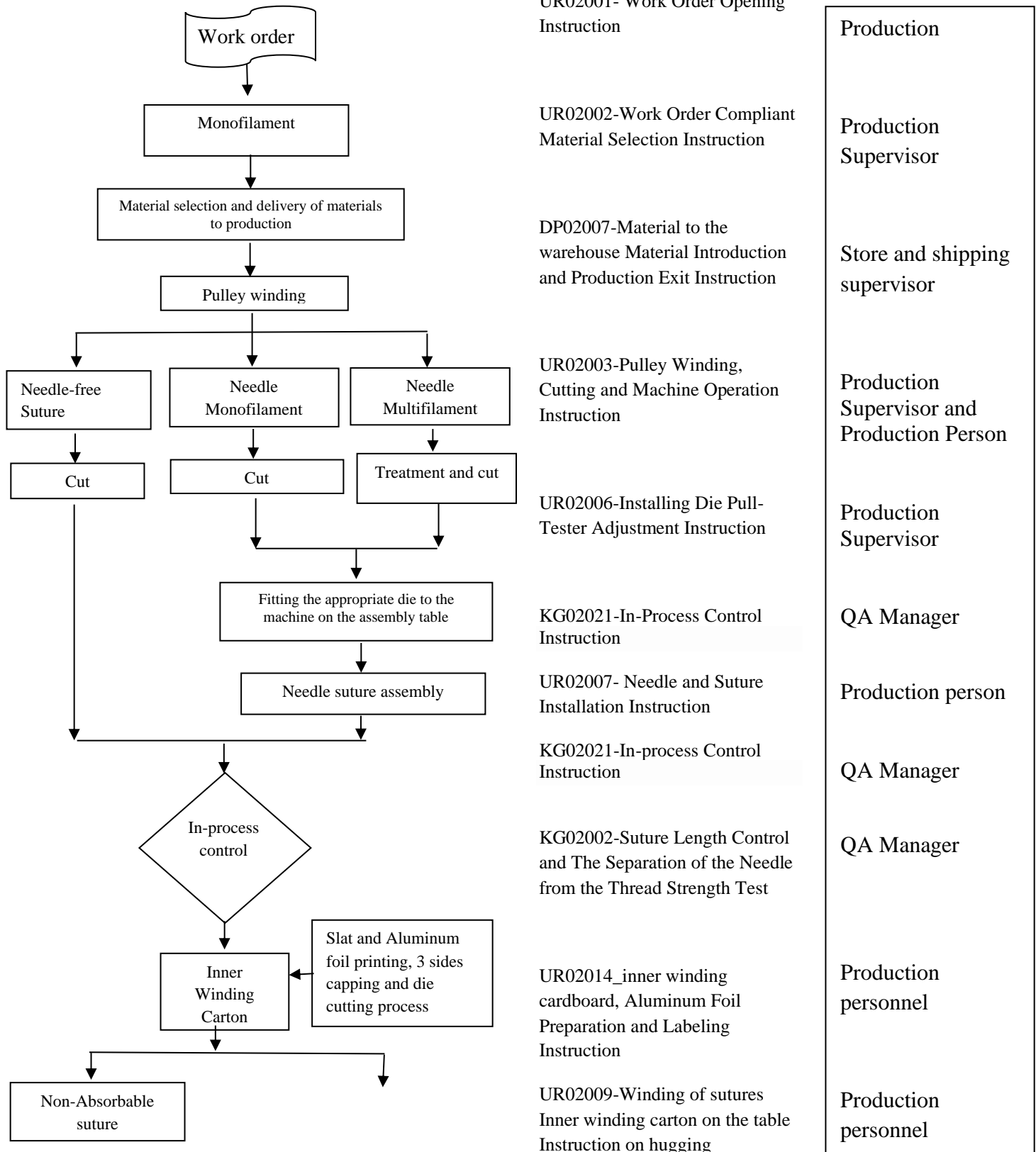

	GMD Group Medikal San. Ve Tic. A.Ş		Ref No.	DC-06
	Department	Device master file	Rev No.	03
	Damacryl and PGA	PGA (Absorbable)	Issue Date	21.07.2023

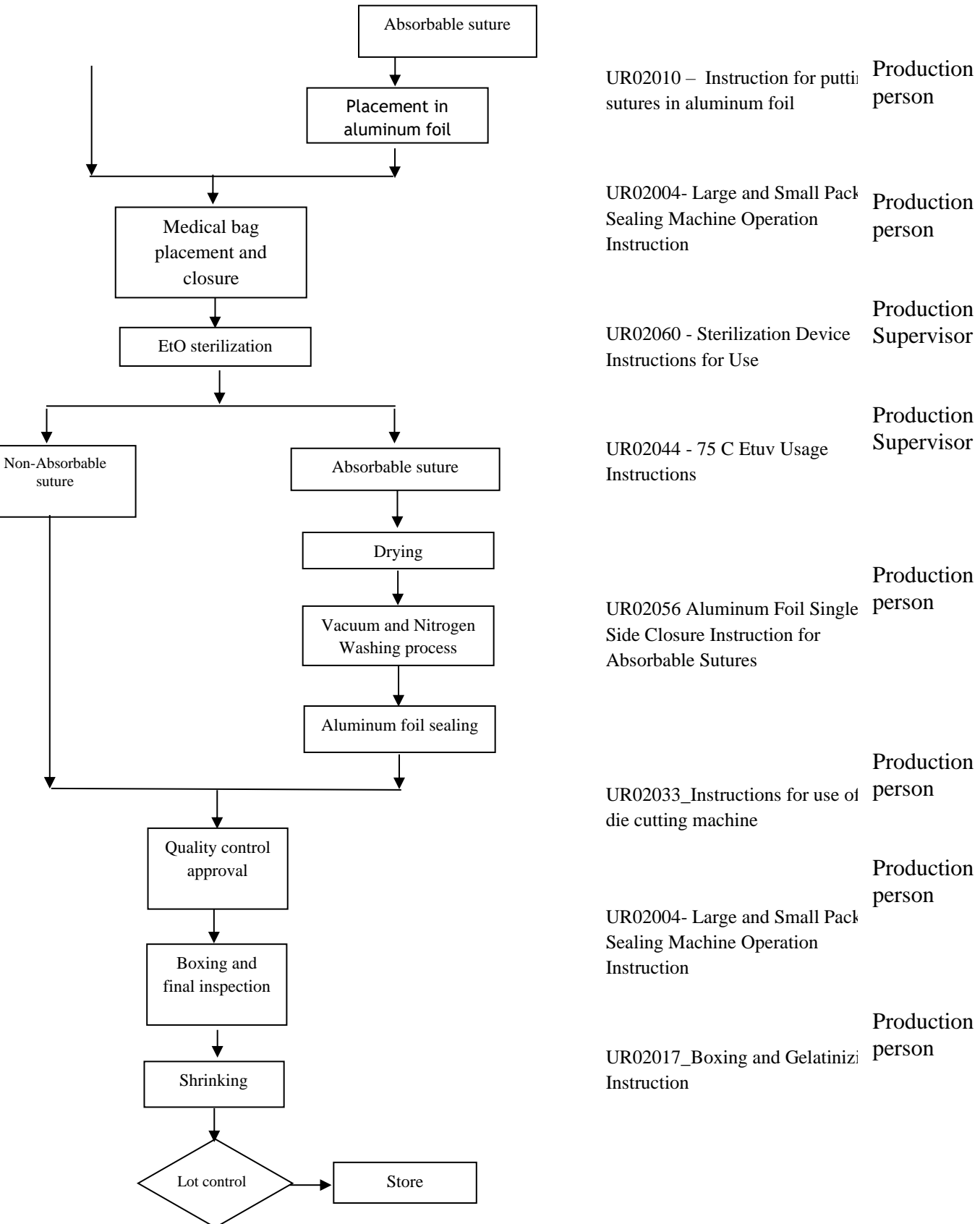

	GMD Group Medikal San. Ve Tic. A.Ş		Ref No.	DC-06
	Department	Device master file	Rev No.	03
	Damacryl and PGA	PGA (Absorbable)	Issue Date	21.07.2023

From material entrance into the production to Temporary warehouse (carrantina area) all processes are carried out in GMD Clean room.

1.19. Product and Test Specifications:

The Damacryl and PGA® range of sutures started to successfully produced and extensively used in 2016(June)

The following tests are specified as final batch release controls according to USP 0: Table:


Test specifications and Lot results

Pilot Production

Test Parameter	Method	Lot 23E02C035 DC0K40	
		Limit	Result
Needle Attachment	USP 0	15,3	34,7
Sterility Test	USP 0	Sterile	Sterile

Validation1

Test Parameter	Method	Lot 16C01 DC1M30	
		Limit	Result
Needle Attachment	USP 1	1,80	5,613
Sterility Test	USP 1	Sterile	Sterile
Water content(moisture)		0.25	0.15

	GMD Group Medikal San. Ve Tic. A.Ş		Ref No.	DC-06
	Department	Device master file	Rev No.	03
	Damacryl and PGA	PGA (Absorbable)	Issue Date	21.07.2023

Validation 2

Parameter	Method	Lot 16C01 DCOM36	
		Limit	Result
Needle Attachment	USP 0	1,80	6,827
Sterility Test	USP 0	Sterile	Sterile
Water content(moisture)		0.25	0.13

Sterility Test	USP 0	Sterile	Sterile	Sterile	Sterile	Sterile	Sterile
Water content(moisture)		0.25	0.15	0.25	0.15	0.25	0.15

The test reports of PGA lots given on the above table are given on TF11 Test reports.

1.20. Stability/Aging Test and Product Life Time :


The relevant tests and reports are given on TF11 Test Reports

1.21. Product Tests(Performance) :

The relevant tests and reports are given on TF11 Test Reports

1.22. The Product Classification

According to annex IX, rule no. 8 of EC-Directive 93/42/EEC the product category Damacryl and PGA is classified as Medical Devices class III. The product is a synthetic absorbable medical device (surgical suture) for surgically invasive long-term implantation applications.

	GMD Group Medikal San. Ve Tic. A.Ş		Ref No.	DC-06
	Department	Device master file	Rev No.	03
	Damacryl and PGA	PGA (Absorbable)	Issue Date	21.07.2023


Rule 8:

“All implantable devices and long-term surgically invasive devices are in Class IIb unless they are intended:

- to be placed in the teeth, in which case they are in Class IIa,
- to be used in direct contact with the heart, the central circulatory system or the central nervous system, in which case they are in Class III,
- to have a biological effect or to be wholly or mainly absorbed, in which case they are in Class III,
- or to undergo chemical change in the body, except if the devices are placed in the teeth, or to administer medicines, in which case they are in Class III.


Biocompatibility Test Matrix

DEVICE CATEGORIES		BIOLOGICAL EFFECT												
		Cytotoxicity	Sensitization	Irritation/Intracutaneous	Acute	Systemic Toxicity	Subchronic Toxicity	Genotoxicity	Implantation	Hemocompatibility	Chronic Toxicity	Carcinogenicity	Reproductive/Developmental	
SURFACE DEVICES	Skin	A	x	x	x									
		B	x	x	x									
		C	x	x	x									
	Mucosal Membrane	A	x	x	x									
		B	x	x	x									
		C	x	x	x		x	x						
	Breached or Compromised Surfaces	A	x	x	x									
		B	x	x	x									
		C	x	x	x		x	x						
Blood Path, Indirect	A	x	x	x	x				x					
	B	x	x	x	x				x					
	C	x	x	x	x	x	x		x					

	GMD Group Medikal San. Ve Tic. A.Ş		Ref No.	DC-06
	Department	Device master file	Rev No.	03
	Damacryl and PGA	PGA (Absorbable)	Issue Date	21.07.2023

In compliance with the standard EN ISO 10993-1, Table 1, the biocompatibility classification of our product and the tests required in accordance are determined as below:

Product	Classification	Required Tests
Damacryl	Implant Devices Tissue/Bone C-Permanent	<ul style="list-style-type: none"> • Cytotoxicity • Sensitization • Irritation • Acute Systemic Toxicity • Subchronic Toxicity • Genotoxicity • Implantation

	GMD Group Medikal San. Ve Tic. A.Ş		Ref No.	DC-06
	Department	Device master file	Rev No.	03
	Damacryl and PGA	PGA (Absorbable)	Issue Date	21.07.2023

The following evaluation is taken by NAMSA, CA. U.S.A.

- Cytotoxicity study using the ISO elution method
 - In vitro hemolysis study
 - ISO sensitization study in the guinea pig
 - ISO acute systemic toxicity study in the mouse
 - ISO acute intracutaneous reactivity study in the rabbit
 - ISO rabbit pyrogen study
 - ISO muscle implantation study in the rabbit (4 weeks)
 - ISO muscle implantation study in the rabbit (26 weeks)
- Genotoxicity: bacterial reverse mutation study.

The results of the biocompatibility tests on the PGA absorbable braid demonstrate the excellent biocompatibility of the PGA material and also demonstrate that the ETO sterilization does not have any negative effects on the polymer material

Based on the clinical experiences with the product PGA, the material documentation as well as the biocompatibility test results on the end product, the material used for the manufacturing of the surgical suture and the final product PGA can be regarded as biocompatible and suitable for a clinical application.

All related test reports, supplier data and related documentation with this section are given on Test Reports of this technical file.

1.23. Biodegradation for Medical Device Damacryl.

Using the polyglycolide and poly(L-lactide) properties as a starting point, it is possible to copolymerize the two monomers to extend the range of homopolymer properties. Copolymers of glycolide with both L-lactide and DL-lactide have been developed for both device and drug delivery applications


Sterility test and absorbibility study in rats by Korea Testing & Research Institute

This study was performed to evaluate sterility and absorbability in rat.


1.24. Labels and Instruction for Use

All information required for use of GMD safely is given on the product label and IFU. The device is intended for professional use only and all the applicable information required according to article 13 of Annex I, 93/42/EEC are given on the label and IFU.

The symbols are selected according to EN ISO 15223:2012 The sample layout of IFU and the product label are given on this technical file.

	GMD Group Medikal San. Ve Tic. A.Ş		Ref No.	DC-06
	Department	Device master file	Rev No.	03
	Damacryl and PGA	PGA (Absorbable)	Issue Date	21.07.2023

Symbols used on the labels:

ÜRÜN ÜZERİNDEKİ SEMBOLLER SYMBOLS USED ON THE PRODUCT	SYMBOLY POUŽITÉ NA PRODUKTU ИСПОЛЬЗУЕМЫЕ ОБОЗНАЧЕНИЯ الرموز المستخدمة على المنتج
 <p>Üretim Tarihi Date Of Manufacture Datum Výroby Срок Изготовления تاريخ التصنيع</p>	 <p>Paket Açılmıy yada Zarar Görmüştü Kullanmayınız Do Not Use If Pack Is Opened or Damaged Nepoužívejte, Pokud Je Balení Otvřené Nebo Poškozené Не Применять При Повреждении Упаковки لا تستخدمها إذا كان الغلاف مفتوح أو متضرر</p>
 <p>Son Kullanma Tarihi Date Of Expiry Datum Spotřebny Использовать До تاريخ نهاية الصلاحية</p>	 <p>Tekrar Kullanmayınız Do-Not Reuse Nepoužívejte Opakovaně Не Использовать Повторно لا تعيد استخدامها</p>
 <p>Sıcaklık Limiti Storage Temperature Limits Teplotní Omezení Ограничение Температуры حدود درجة حرارة التخزين</p>	 <p>Etilen Oksit ile Sterilize Edilmiştir Ethylene Oxide Sterilized Product Produkt Sterilizovaný Ethylenoxidem Простерилизовано Оксидом Этилена منتج معقم بأكسيد الإيثيلين</p>
 <p>Lot Numarası Batch Number Číslo Šarže Код Партии رقم الدفعة</p>	 <p>Tekrar Sterilize Edilmez Do-Not Re-Sterilize Nesterilizujte Znovu Повторно Не Стерилизовать لا تقم بإعادة التعقيم</p>
 <p>Nemden Uzak Tutunuz Avoid Moisture Vyvarujte Se Vlhkosti Бережь От Влага تجنب الرطوبة</p>	 <p>Kullanmadan Önce Kullanım Kılavuzunu Okuyunuz Consult Instructions For Use Viz Návod K Použití Обратитесь К Инструкции По Применению راجع التعليمات للاستخدام</p>
 <p>Güneş Işığına Maruz Bırakmayınız Avoid Direct Sunlight Chraňte Před Přímým Slunečním Zářením Избегайте Прямых Солнечных Лучей تجنب أشعة الشمس المباشرة</p>	 <p>Üretici Firma Manufacturer Výrobce Производитель الشركة المصنعة</p>

1.25. Planned Design Change


Any design change is not planned. If a change is required, it should be carried out in accordance with GE01021 Change Management Procedure

1.26. Sale Information for DAMACRYL and PGA


The product has been in market since 2016. CE marked has been provided from UDEM (2292)

The product has been in the Turkey and European Markets. For 2022 sale informations has been presented below;


DAAMCRYL : PGA 3-0 75CM NEEDLE 60MM STRIHGT CUTTING	BOX	646
--	------------	------------

	GMD Group Medikal San. Ve Tic. A.Ş		Ref No.	DC-06
	Department	Device master file	Rev No.	03
	Damacryl and PGA	PGA (Absorbable)	Issue Date	21.07.2023

DAMACRYL : PGA 0 75CM 45MM CUTTING	BOX	10
DAMACRYL : PGA 0 75CM NEEDLE 48MM ROUND	BOX	100
DAMACRYL : PGA 0 90CM NEEDLE 30MM CUTTING	BOX	24
DAMACRYL : PGA 0 90CM NEEDLE 36MM ROUND	BOX	90
DAMACRYL : PGA 0 Violet 150cm without Needle	BOX	49
DAMACRYL : PGA 0 Violet 75cm Needle Rev.C 3/8 40mm	BOX	599
DAMACRYL : PGA 0 Violet 75cm Needle T.P. 1/2 40mm	BOX	6338
DAMACRYL : PGA 0 Violet 75cm Needle T.P. 1/2 50mm	BOX	524
DAMACRYL : PGA 0 Violet 90cm Needle T.P. 1/2 40mm	BOX	245
DAMACRYL : PGA 0 Violet 90cm Needle T.P. 1/2 50mm	BOX	150
DAMACRYL : PGA 1 100CM NEEDLE 50MM ROUND	BOX	723
DAMACRYL : PGA 1 150CM X 2 without needle	BOX	27
DAMACRYL : PGA 1 150CM X1	BOX	28
DAMACRYL : PGA 1 50CM X7 without needle	BOX	202
DAMACRYL : PGA 1 90CM NEEDLE 30MM ROUND	BOX	53
DAMACRYL : PGA 1 90CM NEEDLE 35MM ROUND	BOX	128
DAMACRYL : PGA 1 Violet 75cm Needle Rev.C 3/8 40mm	BOX	742
DAMACRYL : PGA 1 Violet 75cm Needle T.P. 1/2 40mm	BOX	12140
DAMACRYL : PGA 1 Violet 75cm Needle T.P. 1/2 45mm	BOX	153
DAMACRYL : PGA 1 Violet 75cm Needle T.P. 1/2 50mm	BOX	1351
DAMACRYL : PGA 1 Violet 90cm Needle T.P. 1/2 40mm	BOX	1123
DAMACRYL : PGA 1 Violet 90cm Needle T.P. 1/2 50mm	BOX	363
DAMACRYL : PGA 2 150CM X1	BOX	26
DAMACRYL : PGA 2 150CM X2 without needle	BOX	27
DAMACRYL : PGA 2 50CM X7 WITHOUT NEEDLE	BOX	210
DAMACRYL : PGA 2 Violet 75cm Needle Rev.C 3/8 40mm	BOX	516
DAMACRYL : PGA 2 Violet 75cm Needle T.P. 1/2 40mm	BOX	9780
DAMACRYL : PGA 2 Violet 75cm Needle T.P. 1/2 45mm	BOX	150
DAMACRYL : PGA 2 Violet 75cm Needle T.P. 1/2 50mm	BOX	1386
DAMACRYL : PGA 2 Violet 90cm Needle T.P. 1/2 40mm	BOX	378
DAMACRYL : PGA 2 Violet 90cm Needle T.P. 1/2 50mm	BOX	522
DAMACRYL : PGA 2-0 75CM X 5	BOX	42
DAMACRYL : PGA 2-0 90CM NEEDLE 30MM CUTTING	BOX	26
DAMACRYL : PGA 2-0 90CM NEEDLE 30MM ROUND	BOX	148
DAMACRYL : PGA 2-0 90CM NEEDLE 35MM ROND	BOX	176
DAMACRYL : PGA 2-0 Violet 150cm without Needle	BOX	42
DAMACRYL : PGA 2-0 Violet 75cm Needle Rev.C 3/8 26mm	BOX	864
DAMACRYL : PGA 2-0 Violet 75cm Needle Rev.C 3/8 30mm	BOX	337
DAMACRYL : PGA 2-0 Violet 75cm Needle T.P. 1/2 26mm	BOX	4919
DAMACRYL : PGA 2-0 Violet 75cm Needle T.P. 1/2 30mm	BOX	304

	GMD Group Medikal San. Ve Tic. A.Ş		Ref No.	DC-06
	Department	Device master file	Rev No.	03
	Damacryl and PGA	PGA (Absorbable)	Issue Date	21.07.2023

DAMACRYL : PGA 2-0 Violet 75cm Needle T.P. 1/2 35mm	BOX	330
DAMACRYL : PGA 2-0 Violet 75cm Needle T.P. 1/2 40mm	BOX	1157
DAMACRYL : PGA 3-0 Violet 2x75cm without Needle	BOX	17
DAMACRYL : PGA 3-0 Violet 75cm Needle Rev.C 3/8 26mm	BOX	1583
DAMACRYL : PGA 3-0 Violet 75cm Needle T.P. 1/2 20mm	BOX	207
DAMACRYL : PGA 3-0 Violet 75cm Needle T.P. 1/2 26mm	BOX	2273
DAMACRYL : PGA 3-0 Violet 75cm Needle T.P. 1/2 30mm	BOX	500
DAMACRYL : PGA 4-0 75CM NEEDLE 20MM CUTTING	BOX	236
DAMACRYL : PGA 4-0 75CM NEEDLE 26MM CUTTING	BOX	136
DAMACRYL : PGA 4-0 75CM NEEDLE 26MM ROUND	BOX	2609
DAMACRYL : PGA 4-0 Violet 75cm Needle Rev.C 3/8 16mm	BOX	355
DAMACRYL : PGA 4-0 Violet 75cm Needle T.P. 1/2 16mm	BOX	849
DAMACRYL : PGA 4-0 Violet 75cm Needle T.P. 1/2 20mm	BOX	637
DAMACRYL : PGA 5-0 Violet 75cm Needle Rev.C 3/8 16mm	BOX	96
DAMACRYL : PGA 5-0 Violet 75cm Needle T.P. 1/2 13mm	BOX	220
DAMACRYL : PGA 5-0 Violet 75cm Needle T.P. 1/2 16mm	BOX	132
DAMACRYL : PGA 6-0 Violet 75cm Needle Rev.C 3/8 13mm	BOX	90
DAMACRYL : PGA 6-0 Violet 75cm Needle T.P. 1/2 13mm	BOX	207
DAMACRYL : pga 0 50cm x7 without needle	BOX	67
DAMACRYL :PGA 0 90CM NEEDLE 30MM ROUND	BOX	88
DAMACRYL :PGA 1 75CM X 5	BOX	40
DAMACRYL: PGA 0 100CM NEEDLE 40MM ROUND	BOX	81
DAMACRYL: PGA 0 100CM NEEDLE 50MM ROUND	BOX	79
DAMACRYL: PGA 1 100CM NEEDLE 45MM ROUND	BOX	80
DAMACRYL: PGA 1 100CM NEEDLE 48MM ROUND	BOX	1132
DAMACRYL: PGA 1 150CM Double NEEDLE ROUND 48MM	BOX	43
DAMACRYL: PGA 1 90CM NEEDLE CUTTING 40MM	BOX	99
DAMACRYL: PGA 2 75CM 45CM CUTTING	BOX	50
DAMACRYL: pga 0 90cm needle round 35mm	BOX	347
DAMACRYL: pga 1 90cm needle round 48mm	BOX	175
DAMACRYL: pga 2 90cm 48mm round	BOX	299
DAMACRYL: pga 2-0 90cm needle round 40mm	BOX	265
DAMACRYL:PGA 1 100CM NEEDLE 40MM ROUND	BOX	899
Damacryl : pga 1 75cm needle 36mm cutting	BOX	105
SJ . PGA : PGA 1 75CM NEEDLE 45MM ROUND	BOX	208
SJ PGA : PGA 0 75CM NEEDLE 40MM ROUND	BOX	211
SJ PGA : PGA 0 75CM NEEDLE 45MM ROUND	BOX	207
SJ PGA : PGA 1 75CM NEEDLE 40MM ROUND	BOX	213

	GMD Group Medikal San. Ve Tic. A.Ş		Ref No.	DC-06
	Department	Device master file	Rev No.	03
	Damacryl and PGA	PGA (Absorbable)	Issue Date	21.07.2023

SJ PGA : PGA 2 75CM NEEDLE 40MM ROUND	BOX	214
SJ PGA : PGA 2-0 75CM NEEDLE 26MM ROUND	BOX	93
SJ PGA : PGA 3-0 75CM NEEDLE 75CM 20MM ROUND	BOX	98
SJ PGA : PGA 4-0 75CM NEEDLE 16MM ROUND	BOX	102
TOPLAM		62490

Rev. History

Rev.0 08.07.2016	Initial publish
Rev.1 17.08.2018	Added Product Model Names
Rev.2 25.10.2018	İtems: Intended use of the product,1Intended user, Patient Population, Design Change, Sale Information have been changed and mentioned inthe document by highlight.
Rev.3 12.06.2019	Title update and updated some information

Quality Assurance Manager

Aylin KAYA

