


Lot 7

Iris Retractors

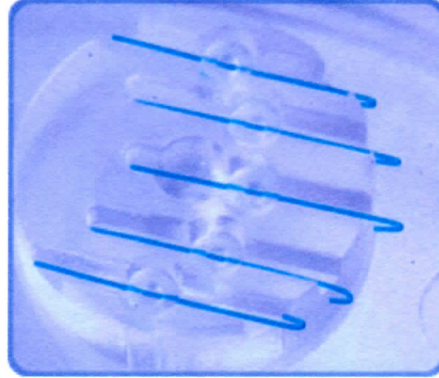

Product's USP:

- Available as a pack of 5
- Soft contoured for comfort
- ETO Sterilization

Best Suitable: For Patients with non-dilating iris conditions.

Indications for Use :

Also known as Iris hooks, these are used in patients having a non-dilating or iris syndromes or in conditions where dilation thru other means like pharma applications is unsuccessful. It is used for mechanical dilation of the pupil and helps to stretch the iris during a surgery and thereby reduce the surgical time.

Another Quality Product from :



Action Medical Mktg Pvt Ltd.

(An ISO 9001:2008 & ISO 13485:2003 Certified Company)

Gat No 528, Koregaon-Bhima, Near SBI, Pune-Nagar Road

Taluka Shirur, Dist Pune Pin 412216, INDIA

Tel : 0091-2137-253617 / 18. Fax : 0091-2137-253620.

E-mail : action@actionmedihealth.com



Website : www.actionmedihealth.com

CATARACT & IOL

IrisCare

Flexible Iris Retractors

Flexible iris retractors are used for temporary retraction (pupil enlargement) of the iris where dilation cannot be achieved pharmacologically.

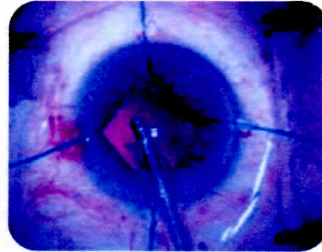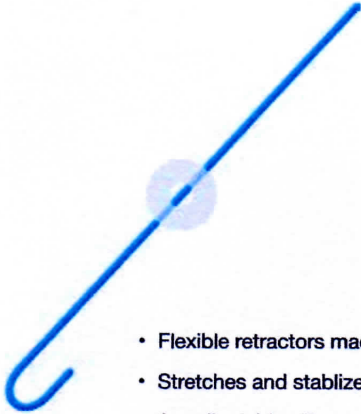

- Flexible retractors made of bright blue PP (Polypropylene)
- Stretches and stabilizes the Iris during the cataract surgery
- An adjustable silicon stopper provides a superior grip and allows adjustment of the pupil size
- Easy to insert and remove
- CE certified


Product Code	Description	Qty. per box
MIPL/D1	Iris Retractors	5 pcs. (Sterile)

Lot 7

CapsuleCare

A Capsular Bag Supporting Device

Capsule hooks are used for temporary support of the capsular bag during a cataract surgery when the zonules are weak or in case of subluxated lenses.

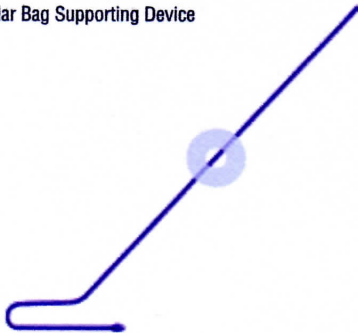

Lot 6

- Capsule hooks are made of Nylon with an adjustable silicon lock
- Designed such that the hook is angled to be in the plane of the anterior capsule
- Easy to get a firm and reliable grip on the rim
- Tip of the hook is rounded to prevent any accidental damage to the peripheral capsule while placing the hooks in the bag
- Round tip also provides a larger support area per hook there by reducing the number of hooks required in each surgery
- Silicon lock helps in retaining the hook in retracted position during the surgery
- CE certified

Product Code	Description	Qty. per box
MIPL/D8	Capsule Hooks	5 pcs. (Sterile)

Advantages

- Maintains the position of the capsule during phaco
- Supports the fornix as well as the capsulorhexis rim, so that the peripheral capsule isn't attracted towards the phaco tip
- Avoids the problem of the capsular tension ring trapping cortex
- Prevents zonular stress during nuclear rotation
- Makes it possible to perform posterior chamber phaco even with subluxated lenses having minimum zonular support
- Easy to insert and remove