

RELIABLE HEARING IN ANY SITUATION

PRODUCT CATALOG



REXTON


OUR LATEST INNOVATIONS FOR RELIABLE HEARING



MOTION CORE TECHNOLOGY



MOTION CORE HEARING AIDS

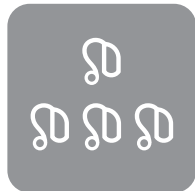

HEARING AID PORTFOLIO



M-CORE R

Updated design, with Rechargeable, CROS, and Telecoil options

6



M-CORE SR

Award-winning design – now with exchangeable receivers

7



M-CORE IX

Our smallest instant-fit CIC, now with Motion Core

8

MYCORE HEARING AIDS

10

TRUCORE HEARING AIDS

40

ACCESSORIES

60

MOTION CORE TECHNOLOGY



Today's hearing aids analyze hearing situations solely on acoustic criteria, then classify them into a small range of predefined "acoustic environments". But when people rely on their senses to keep them safe and sound, hearing clearly in every situation is a must – especially when they're on the job and moving through the day, where environments can change quickly.

Our latest range of hearing aids use proven Motion Core technology to automatically recognize and adapt to more hearing situations than ever. This ensures users automatically hear speech clearer and get the most reliable impression of their environment. M-Core hearing aids are designed to be durable, practical and easy-to-use, so your customers can rely on Rexton in any situation.

UP TO 20 TIMES MORE HEARING SITUATIONS THAN BEFORE, AUTOMATICALLY



HIGH RESOLUTION SOUNDPRO

Using a combination of inputs including real-time motion data, signal to noise ratio, and sound directionality, High Resolution Soundpro can continuously adjust the sound settings to provide the most reliable impression of the user's environment, automatically.



MOTION SENSOR

Motion Core technology features an integrated real-time motion sensor that detects a user's movement, and uses the data to even more accurately determine their sound situation. With this, M-Core hearing aids can automatically detect up to 127 different sound situations.



MY VOICE

The hearing aid recognizes the user's specific voice pattern, and integrates it naturally into what they hear. In turn, this provides the user with a natural-sounding voice that leads to high acceptance rates. And for that the user will communicate more and be more socially active, greatly improving their quality of life.



STREAMING

Audio signals from smartphones, TVs or devices can be streamed wirelessly to the hearing aids; directly or via transmitter. Direct Audio Streaming is available for iOS smartphones, and Android smartphones can stream with the addition of the Smart Mic accessory.



A HEARING AID WEARER'S ALL-IN-ONE SMARTPHONE COMPANION

The Rexton App allows wearers to control their hearing aids conveniently and discreetly, e.g. change programs, volume, tone.

Download the Rexton App now!





A FULL RANGE OF RELIABLE HEARING OPTIONS FOR YOUR CUSTOMERS


	Instant-fit					Custom				Slim-RIC				RIC					BTE				
	M-CORE iX-CIC (inoX)	Sterling IIC	Sterling CIC	Sterling ITC	Sterling ITE	M-CORE SR	StyleLine	StyleLine BT	Emerald XS	M-Core R312 (Emerald S)	M-Core R-Li (Stellar RIC)	Emerald M	Stellar	Mosaic S	Mosaic M	Stellar M	Mosaic P	Mosaic HP					
Battery size	10 / 312	10	10	312	312	Rechargeable battery	Rechargeable battery	Rechargeable battery	10	312	Rechargeable battery	13	Rechargeable battery	312	13	Rechargeable battery	13	675					
PL 80	<div><div></div><div></div><div>New</div></div>	<div><div></div></div>	<div><div></div><div></div></div>	<div><div></div><div>New</div></div>	<div><div></div><div>New</div></div>	<div><div></div><div>New</div></div>	<div><div></div></div>	<div><div></div></div>	<div><div></div></div>	<div><div></div><div></div></div>	<div><div></div><div>New</div></div>	<div><div></div><div>New</div></div>	<div><div></div><div></div></div>	<div><div></div></div>		<div><div></div><div></div></div>	<div><div></div></div>	<div><div></div><div></div></div>	<div><div></div></div>				
PL 60	<div><div></div><div></div><div>New</div></div>	<div><div></div></div>	<div><div></div><div></div></div>	<div><div></div><div>New</div></div>	<div><div></div><div>New</div></div>	<div><div></div><div>New</div></div>	<div><div></div></div>	<div><div></div></div>	<div><div></div><div></div></div>	<div><div></div><div>New</div></div>	<div><div></div><div>New</div></div>	<div><div></div><div></div></div>	<div><div></div></div>		<div><div></div><div></div></div>	<div><div></div></div>	<div><div></div><div></div></div>	<div><div></div></div>					
PL 40	<div><div></div><div></div><div>New</div></div>	<div><div></div></div>	<div><div></div><div></div></div>	<div><div></div><div>New</div></div>	<div><div></div><div>New</div></div>	<div><div></div><div>New</div></div>	<div><div></div></div>	<div><div></div></div>	<div><div></div><div></div></div>	<div><div></div><div>New</div></div>	<div><div></div><div>New</div></div>	<div><div></div><div></div></div>	<div><div></div></div>	<div><div></div></div>	<div><div></div><div></div></div>	<div><div></div></div>	<div><div></div><div></div></div>	<div><div></div></div>					
PL 30	<div><div></div><div></div><div>New</div></div>		<div><div></div></div>	<div><div></div></div>	<div><div></div></div>		<div><div></div></div>		<div><div></div><div>New</div></div>	<div><div></div><div></div></div>	<div><div></div></div>	<div><div>New</div></div>	<div><div></div></div>	<div><div></div></div>	<div><div></div><div>New</div></div>	<div><div></div></div>	<div><div></div><div>New</div></div>	<div><div></div></div>					
PL 20	<div><div></div><div>New</div></div>		<div><div></div></div>	<div><div></div></div>	<div><div></div></div>		<div><div></div></div>		<div><div></div><div>New</div></div>	<div><div></div><div></div></div>	<div><div></div></div>	<div><div>New</div></div>		<div><div></div></div>	<div><div></div><div>New</div></div>	<div><div></div></div>	<div><div></div><div>New</div></div>	<div><div></div></div>					
PL 10	<div><div></div></div>		<div><div></div></div>	<div><div></div></div>	<div><div></div></div>					<div><div></div></div>				<div><div></div></div>	<div><div></div></div>		<div><div></div></div>	<div><div></div></div>					
	<div><div></div><div></div></div>	<div><div></div><div></div></div>	<div><div>New</div></div>	<div><div></div></div>	<div><div>New</div></div>																		

OUR MOST TRUSTED RIC

NOW WITH MOTION CORE TECHNOLOGY

Our most trusted RIC hearing aid now comes with proven Motion Core technology that automatically recognizes and adapts to more hearing situations than ever. This ensures users automatically hear speech clearer and get the most reliable impression of their environment. M-Core RIC hearing aids are designed to be durable, practical and easy-to-use – with either 312 battery or rechargeable option.

Available colors

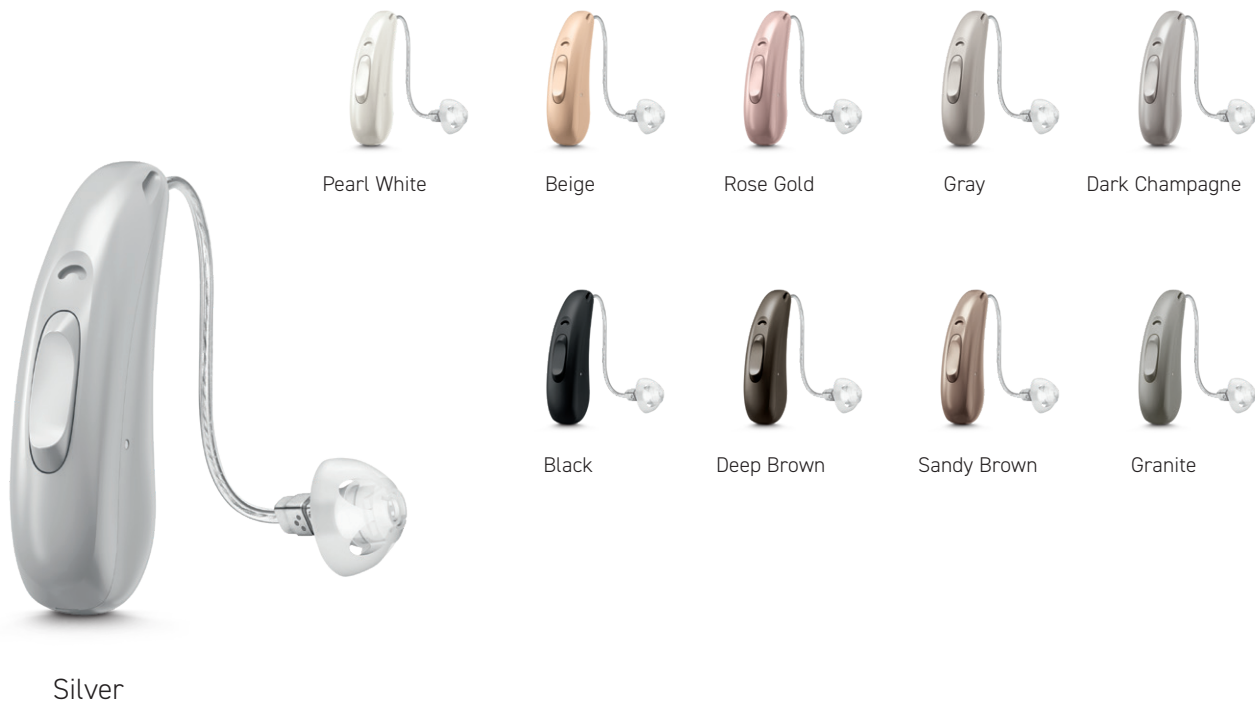

M-CORE R-Li



Performance Level	80, 60, 40
Battery	Lithium-ion battery
Accessories	Charging Station
	Smart Mic
	Smart Transmitter 2,4
	Rexton App
Features	Integrated Motion Sensor
	Programmable rocker switch
	Mini Receivers 2.0
	Color Conversion Kits
	Programming via: > NoahLink Wireless

- M-Core R-Li features a rechargeable lithium-ion battery format. The rechargeable batteries have a running time of up to 21 h incl. 5 h streaming on one full charge, and the quick charge feature delivers up to 6 hours extra running time in 30 minutes



Connexx Charging Case

M-CORE R312



Performance Level	80, 60, 40
Battery	312
Accessories	Smart Mic
	Smart Transmitter 2,4
	Rexton App
Features	Integrated Motion Sensor
	Programmable rocker switch
	Mini Receivers 2.0
	Color Conversion Kits
	Color Conversion Kits with design integrated telephone coil (M-Core R312T)
	Programming via: > NoahLink Wireless

- M-Core R312 is optionally available with an integrated Telecoil for reliable hearing where loop systems are available



M-Core R-312T with optional telecoil



FITTING RANGES

MOTION CORE RIC HEARING











M-CORE R-LI



M-CORE R312



Battery size	Rechargeable battery (Li-Ion)	312
--------------	-------------------------------	-----

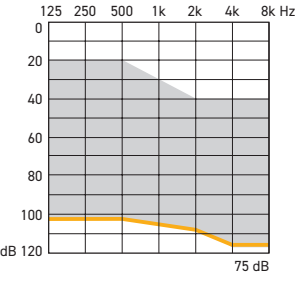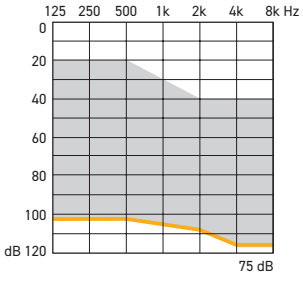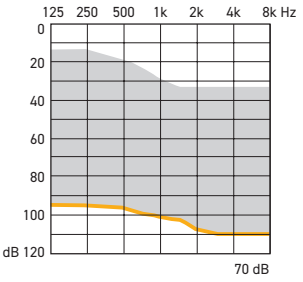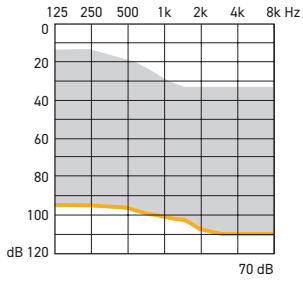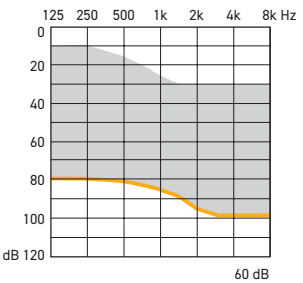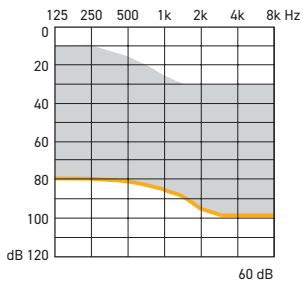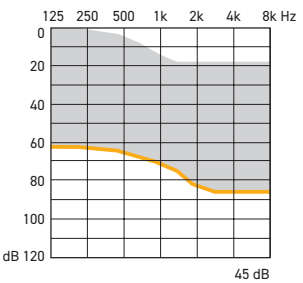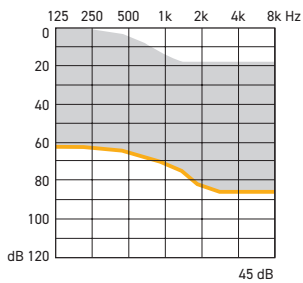

M-CORE R312



Performance Level 80 | 60 | 40

New



Highlights

- › Integrated motion sensor
- › Battery size 312
- › Programmable rocker switch
- › Mini Receivers 2.0
- › Smart Demo device available
- › Color Conversion Kits with completely integrated telephone coil
- › Programming via: NoahLink Wireless

THE ENTIRELY NEW DESIGNED UNIVERSAL DEVICE

Taking hearing to a new level. With the new MOTION CORE technology, we set a new standard for hearing comfort. With over 120 hearing situations M-Core R312 automatically offers the right setting in every situation, thanks to the integrated motion sensor. Optionally available with a completely integrated telephone coil which can be easily fitted with a color conversion kit.

Receiver	S	M	P	HP
FOG/OSPL in dB	45/108	60/119	70/124	75/130
Battery current drain	1.7 mA	1.9 mA	1.8 mA	1.8 mA

Colors



Features / Performance Level	80	60	40
Channels/Controls/Programs	48/20/6	32/16/6	24/12/6
Soundpro	HighRes	HighRes	HighRes
My Voice	•	•	•
Direct Streaming/ Auto Volume	Made for iPhone, Android via Smart Mic/Auto Volume	Made for iPhone, Android via Smart Mic/Auto Volume	Made for iPhone, Android via Smart Mic/Auto Volume
Wireless sync	•	•	•
Directionality	Automatic / Adaptive iOmni Front & Back Left & Right Narrow	Automatic / Adaptive iOmni Front & Back	Automatic / Adaptive iOmni
Noise Reduction	Noise Management SoundSmoothing Directional	Noise Management SoundSmoothing Directional	Noise Management SoundSmoothing Directional
Wind Noise Reduction	Standard Binaural	Standard Binaural	Standard
Reverb Reducer	•	-	-
Bandwidth: Extension / Compression	• / •	- / •	- / •
Music Enhancer: Live / Recording / Acting	•	•	-
Tinnitus	Sound Therapy Notch Therapy	Sound Therapy Notch Therapy	Sound Therapy Notch Therapy
Xphone	•	•	•
Acclimatization/Data Logging	• / •	• / •	• / •
Accessories & Apps			
Smart Mic	•	•	•
Smart Transmitter 2,4	•	•	•
Smart Key	•	•	•
Rexton App	•	•	•

M-CORE R-Li



Performance Level 80 | 60 | 40

New



Highlights

- › Integrated motion sensor
- › Lithium-ion battery
- › Programmable rocker switch
- › Mini Receivers 2.0
- › Smart Demo device available
- › Color Conversion Kits
- › Programming via: NoahLink Wireless

THE NEW ULTIMATE RIC DEVICE

Our best sound impression ever! With M-Core R-Li we combine the best hearing comfort, thanks to over 120 hearing situations with a new housing which is more than 15% smaller than its predecessor. All powered with a strong, environment friendly lithium-ion battery with up to 23h runtime. Of course with all connectivity options you value from our Rexton portfolio.

Receiver	S	M	P	HP
FOG/OSPL in dB	45/108	60/119	70/124	75/130
Batt. runtime w/o streaming	up to 23 h	up to 23 h	up to 23 h	up to 23 h
Batt. runtime w. 5h streaming	up to 21 h	up to 20 h	up to 20 h	up to 20 h

Colors



Features / Performance Level	80	60	40
Channels/Controls/Programs	48/20/6	32/16/6	24/12/6
Soundpro	HighRes	HighRes	HighRes
My Voice	•	•	•
Direct Streaming/ Auto Volume	Made for iPhone, Android via Smart Mic/Auto Volume	Made for iPhone, Android via Smart Mic/Auto Volume	Made for iPhone, Android via Smart Mic/Auto Volume
Wireless sync	•	•	•
Directionality	Automatic / Adaptive iOmni Front & Back Left & Right Narrow	Automatic / Adaptive iOmni Front & Back Narrow	Automatic / Adaptive iOmni Narrow
Noise Reduction	Noise Management SoundSmoothing Directional	Noise Management SoundSmoothing Directional	Noise Management SoundSmoothing Directional
Wind Noise Reduction	Standard Binaural	Standard Binaural	Standard
Reverb Reducer	•	-	-
Bandwidth: Extension / Compression	• / •	- / •	- / •
Music Enhancer: Live / Recording / Acting	•	•	-
Tinnitus	Sound Therapy Notch Therapy	Sound Therapy Notch Therapy	Sound Therapy Notch Therapy
Xphone	•	•	•
Acclimatization/Data Logging	• / •	• / •	• / •
Accessories & Apps			
Connexx charging station	•	•	•
Smart Mic	•	•	•
Smart Transmitter 2,4	•	•	•
Smart Key	•	•	•
Rexton App	•	•	•

AWARD-WINNING DESIGN

NOW WITH EXCHANGEABLE RECEIVERS

Our award-winning Slim-RIC design now features proven Motion Core technology that automatically recognizes and adapts to up to 20 times more hearing situations than before. This ensures users hear speech clearer and get the most reliable impression of their environment, automatically. And now with exchangeable receivers, you can easily fit M-Core SR hearing aids to more users than ever.

Choose from 5 distinct color combinations



M-CORE SR



Performance Level	80, 60, 40
Battery	Lithium-ion battery
Accessories	Portable charger
	Smart Mic
	Smart Transmitter 2,4
	Rexton App
Features	Integrated Motion Sensor
	127 sound situations
	Exchangeable Receivers
	Direct Streaming
	Programming via: > NoahLink Wireless

Updated portable charger for power on-the-go

M-Core SR hearing aids feature rechargeable lithium-ion batteries for hassle-free power. The included charging case enables on-the-go top-ups, for up to 4 days of use – so your customers can enjoy a relaxing weekend without needing a power socket.



M-Core SR

Connexx portable charging case

INSTANT FIT

WITH THE HIDDEN POWER OF MOTION CORE



M-CORE iX

Performance Level	80, 60, 40
Battery	10
Accessories	Smart Key
	Rexton App
Features	Dynamic Soundscape Processing
	70 sound situations
	Instant-Fit Sleeves
	Programming via: Flex connector > HiPro, Connexx Air

ACCESSORIES





FITTING RANGES

SLIM-RIC AND INSTANT FIT CIC HEARING AIDS





M-CORE SR



M-CORE IX



Battery size	Rechargeable battery (Li-Ion)	10
--------------	-------------------------------	----

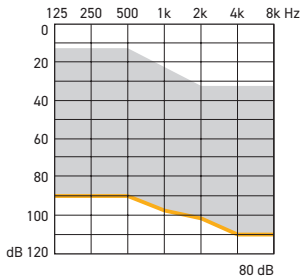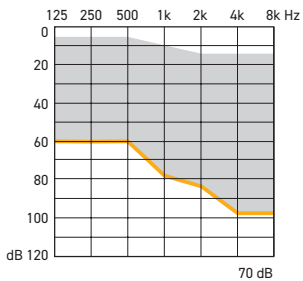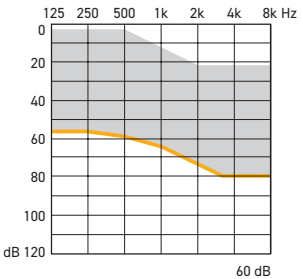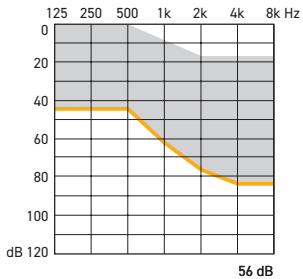

M-CORE SR



Performance Level 80 | 60 | 40



AWARD-WINNING DESIGN
- NOW WITH EXCHANGEABLE
RECEIVERS

Our award-winning Slim-RIC design now features proven Motion Core technology that automatically recognizes and adapts to up to 20 times more hearing situations than before. This ensures users hear speech clearer and get the most reliable impression of their environment, automatically. And now with exchangeable receivers, you can easily fit M-Core SR hearing aids to more users than ever.

Highlights

- › Lithium-ion battery
- › Smart Demo device available
- › Programming via:
NoahLink Wireless

Receiver	S	M	P
FOG/OSPL in dB	45/108	60/119	70/124
Battery runtime	up to 19 h	up to 19 h	up to 19 h

Colors



Features / Performance Level	80	60	40
Channels/Controls/Programs	48/20/6	32/16/6	24/12/6
HD Soundpro	•	•	•
My Voice	•	•	•
Direct Streaming/ Auto Volume	• / •	• / •	• / •
Wireless sync	•	•	•
Directionality	Automatic/Adaptive iOmni Front/Back Left/Right Narrow	Automatic/Adaptive iOmni Front/Back	Automatic/Adaptive iOmni
Noise Reduction	Noise Management SoundSmoothing Directional	Noise Management SoundSmoothing Directional	Noise Management SoundSmoothing
Wind Noise Reduction	Standard Binaural	Standard Binaural	Standard
Reverb Reducer	•	-	-
Bandwidth: Extension / Compression	• / •	- / •	- / •
Music Enhancer: Live / Recording / Acting	• / • / •	• / • / •	- / - / •
Tinnitus	Sound Therapy Notch Therapy	Sound Therapy Notch Therapy	Sound Therapy Notch Therapy
Xphone	•	•	• / - / -
Acclimatization/Data Logging	Premium / •	Premium / •	High / •
Accessories & Apps			
Mobile charging system	•	•	•
Smart Key	•	•	•
Rexton App	•	•	•
Smart Mic	•	•	•
Smart Transmitter 2,4	•	•	•

M-CORE iX



Performance Level 80 | 60 | 40



INSTANT FIT WITH THE HIDDEN
POWER OF MOTION CORE


Our virtually invisible CIC design now features proven Motion Core technology that automatically recognizes and adapts to up to 10 times more hearing situations than before. This ensures users hear speech clearer and get the most reliable impression of their environment, automatically. M-Core iX hearing aids are designed to be discreet, practical and easy-to-use.

Highlights

- › 10 Size battery
- › Smart Demo device available
- › Programming via:
ConnexxAir
Flex connector

FOG/OSPL in dB	60/125
----------------	--------

Colors

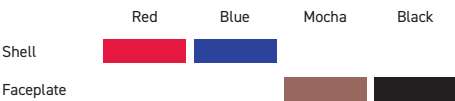

Features / Performance Level	80	60	40
Channels/Controls/Programs	48/20/6	32/16/6	24/12/6
Soundpro	-	-	-
My Voice	-	-	-
Direct Streaming/ Auto Volume	-	-	-
Wireless sync	•	•	•
Directionality	Binaural oneMic Directionality	Binaural oneMic Directionality	Binaural oneMic Directionality
Noise Reduction	Noise Management SoundSmoothing Directional	Noise Management SoundSmoothing Directional	Noise Management SoundSmoothing
Wind Noise Reduction	Standard Binaural	Standard Binaural	Standard
Reverb Reducer	•	-	-
Bandwidth: Extension / Compression	• / •	- / •	- / •
Music Enhancer: Live / Recording / Acting	• / • / •	• / • / •	-
Tinnitus	Sound Therapy Notch Therapy	Sound Therapy Notch Therapy	Sound Therapy Notch Therapy
Xphone	•	•	•
Acclimatization/Data Logging	• / •	• / •	• / •
Accessories & Apps			
Rexton App	•	•	•
Smart Key	•	•	•

MYCORE

A PERFECT COMBINATION OF NATURAL SOUND
AND DIRECT STREAMING



DIRECT STREAMING

Optimally connected to smartphones, TV's and more



MY VOICE

Natural perception of wearer's own voice



STEREO ILOCK

Improves the signal-to-noise ratio in difficult listening situations



VOICE RANGER

Easier one-on-one listening in noisy situations



DYNAMIC EXTENDER

Automatically optimizes the dynamic bandwidth of sounds

FITTING RANGES

CUSTOM HEARING AIDS



INOX



Battery size

10

STERLING IIC & CIC

IIC ONLY AVAILABLE MATRIX
113/50



10

STERLING ITC












312

STERLING ITE



312

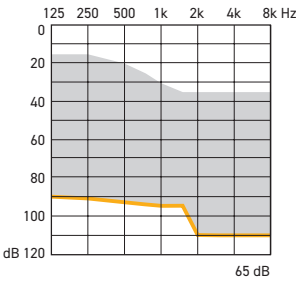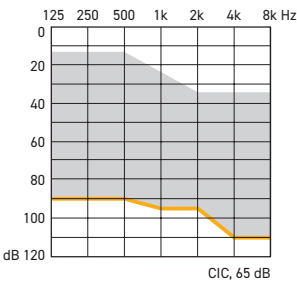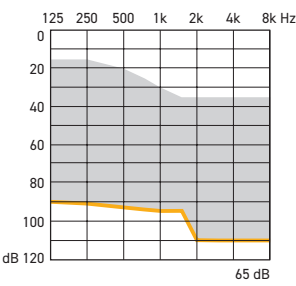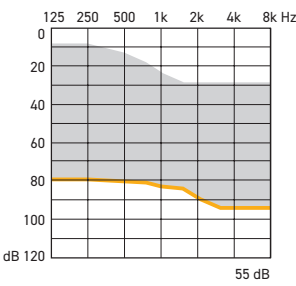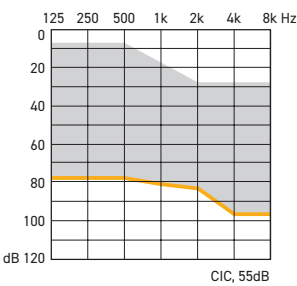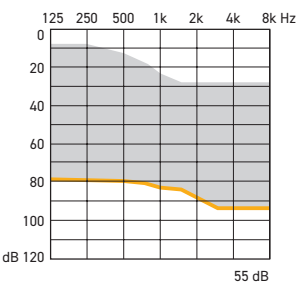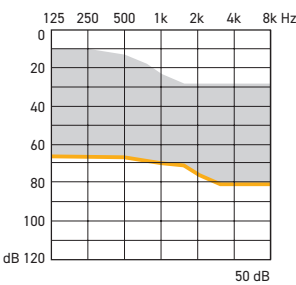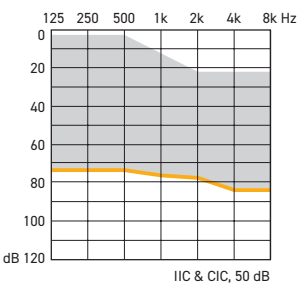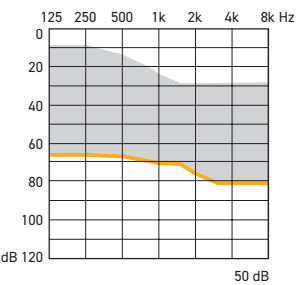



DISCREET AND READY TO FIT

inoX combines the discretion of a custom model without having to take an ear impression. It is not only 20% smaller than its predecessor but also offers enormous sound thanks to the MyCore technology. Even directional characteristics are imbedded in a binaural fitting. With its high flexibility inoX can be fitted to nearly every ear.




Highlights

- › Battery size 10
- › Smart Demo device available
- › Programming via:
- › Flex connector › HiPro
- › ConnexxAir

FOG/OSPL in dB	50/113
Battery current drain	1.1 mA

Colors

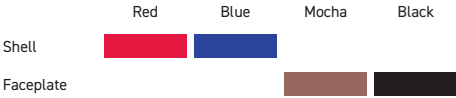

Features / Performance Level	80	60	40	30	20
Channels/Controls/Programs	48/20/6	32/16/6	24/12/6	16/8/4	16/8/4
Soundpro	•	•	•	-	-
My Voice	-	-	-	-	-
Voice Ranger	•	•	•	•	•
Direct Streaming/ Auto Volume	- / -	- / -	- / -	- / -	- / -
Wireless sync	•	•	•	•	•
Directionality	Automatic / Adaptive bin. Single-mic Directionality	Automatic / Adaptive bin. Single-mic Directionality	Automatic / Adaptive bin. Single-mic Directionality	-	-
Noise Reduction	Noise Management SoundSmoothing Directional	Noise Management SoundSmoothing Directional	Noise Management SoundSmoothing	Noise Management SoundSmoothing	Noise Management
Wind Noise Reduction	Standard	Standard	Standard	Standard	-
Reverb Reducer	•	-	-	-	-
Bandwidth: Extension / Compression	• / •	- / •	- / •	- / •	- / •
Music Enhancer: Live / Recording / Acting	• / • / •	• / - / -	-	-	-
Tinnitus	Sound Therapy Notch Therapy	Sound Therapy Notch Therapy	Sound Therapy Notch Therapy	Sound Therapy	-
Xphone	•	•	•	-	-
Acclimatization/Data Logging	Premium / •	Premium / •	High / •	High / •	Standard / •

Accessories & Apps					
Smart Key	•	•	•	•	•
CROS inoX 8C	•	•	•	-	-
Rexton App	•	•	•	•	•



ULTRA-DISCREET
WITH MYCORE FEATURES

Our smallest custom models to offer customers discretion and performance. Maximum discretion ensured by small components and custom-made housing. Binaural signal processing ensures easier listening in loud environments. Discreet Smart key or Rexton App allows wearers to change volume, programs and more with discretion.




Highlights

- › Battery size 10
- › Programming via:
- › ConnexxAir
- › FlexConnector (CIC only)


	50 dB	55db (CIC only)	65 dB (CIC only)
FOG/OSPL in dB	50/113	55/118	65/124
Battery current drain	1.3 mA	1.3 mA	1.3 mA

Colors

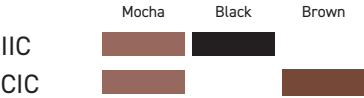

Features / Performance Level	80	60	40
Channels/Controls/Programs	48/20/6	32/16/6	24/12/6
Soundpro	•	•	•
My Voice	-	-	-
Voice Ranger	•	•	•
Direct Streaming/ Auto Volume	- / -	- / -	- / -
Wireless sync	•	•	•
Directionality	Automatic / Adaptive bin. Single-mic Directionality	Automatic / Adaptive bin. Single-mic Directionality	Automatic / Adaptive bin. Single-mic Directionality
Noise Reduction	Noise Management SoundSmoothing Directional	Noise Management SoundSmoothing Directional	Noise Management SoundSmoothing
Wind Noise Reduction	Standard	Standard	Standard
Reverb Reducer	•	-	-
Bandwidth: Extension / Compression	• / •	- / •	- / •
Music Enhancer: Live / Recording / Acting	• / • / •	• / - / -	-
Tinnitus	Sound Therapy Notch Therapy	Sound Therapy Notch Therapy	Sound Therapy Notch Therapy
Xphone	•	•	•
Acclimatization/Data Logging	Premium / •	Premium / •	High / •

Accessories & Apps			
Smart Key	•	•	•
CROS inoX 8C	•	•	•
Rexton App	•	•	•

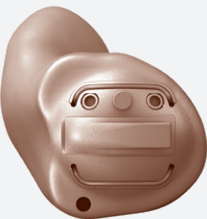

Highlights

- › Battery size 312
- › T-coil option (non-BT)
- › Programming via:

› NoahLink Wireless (BT)

› ConnexxAir / › HiPro

	50 dB	55 dB	65 dB
FOG/OSPL in dB	50/113	55/118	65/124
Battery current drain	1.3 mA	1.3 mA	1.3 mA

Colors

Beige

Brown


Tan

Dark Brown

Mocha (non-BT)

Features / Performance Level	80 <div>New</div>	80 BT	60 <div>New</div>	60 BT	40 <div>New</div>	40 BT
Channels/Controls/Programs	48/20/6	48/20/6	32/16/6	32/16/6	24/12/6	24/12/6
Soundpro	•	•	•	•	•	•
My Voice	-	-	-	-	-	-
Voice Ranger	•	•	•	•	•	•
Direct Streaming/ Auto Volume	-	Made for iPhone via TV Transmitter & Smart Mic / Auto Volume	-	Made for iPhone via TV Transmitter & Smart Mic / Auto Volume	-	Made for iPhone via TV Transmitter & Smart Mic / Auto Volume
Wireless sync	•	•	•	•	•	•
Directionality	Automatic / Adaptive Front & Back Left & Right Narrow - Stereo iLock	Automatic / Adaptive Front & Back Left & Right Narrow - Stereo iLock	Automatic / Adaptive Front & Back Narrow - iLock	Automatic / Adaptive Front & Back Narrow - iLock	Automatic / Adaptive Narrow - iLock	Automatic / Adaptive Narrow - iLock
Noise Reduction	Noise Management SoundSmoothing Directional	Noise Management SoundSmoothing Directional	Noise Management SoundSmoothing Directional	Noise Management SoundSmoothing Directional	Noise Management SoundSmoothing	Noise Management SoundSmoothing
Wind Noise Reduction	Standard	Standard Binaural	Standard	Standard Binaural	Standard	Standard
Reverb Reducer	•	•	-	-	-	-
Bandwidth: Extension / Compression	• / •	• / •	- / •	- / •	- / •	- / •
Music Enhancer: Live / Recording / Acting	• / • / •	• / • / •	• / - / -	• / - / -	-	-
Tinnitus	Sound Therapy Notch Therapy	Sound Therapy Notch Therapy	Sound Therapy Notch Therapy	Sound Therapy Notch Therapy	Sound Therapy Notch Therapy	Sound Therapy Notch Therapy
Xphone	•	•	•	•	•	•
Acclimatization/Data Logging	Premium / •	Premium / •	Premium / •	Premium / •	High / •	High / •

Accessories & Apps		
Smart Mic	-	•
Smart Transmitter 2,4	-	•
Smart Key	•	•
CROS inoX 8C	•	-
CROS RIC 8C	-	•
Rexton App	•	•

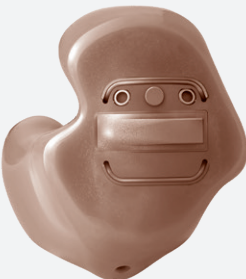

Highlights

- › Battery size 312
- › T-coil option (non-BT)
- › Programming via:

› NoahLink Wireless (BT)

› ConnexxAir / › HiPro

	55 dB	65 dB
FOG/OSPL in dB	55/118	65/124
Battery current drain	1.3 mA	1.3 mA

Colors

Beige

Brown

Tan

Dark Brown

Mocha (non-BT)

Features / Performance Level	80 <div>New</div>	80 BT	60 <div>New</div>	60 BT	40 <div>New</div>	40 BT
Channels/Controls/Programs	48/20/6	48/20/6	32/16/6	32/16/6	24/12/6	24/12/6
Soundpro	•	•	•	•	•	•
My Voice	-	-	-	-	-	-
Voice Ranger	•	•	•	•	•	•
Direct Streaming/ Auto Volume	-	Made for iPhone via TV Transmitter & Smart Mic / Auto Volume	-	Made for iPhone via TV Transmitter & Smart Mic / Auto Volume	-	Made for iPhone via TV Transmitter & Smart Mic / Auto Volume
Wireless sync	•	•	•	•	•	•
Directionality	Automatic / Adaptive Front & Back Left & Right Narrow - Stereo iLock	Automatic / Adaptive Front & Back Left & Right Narrow - Stereo iLock	Automatic / Adaptive Front & Back Narrow - iLock	Automatic / Adaptive Front & Back Narrow - iLock	Automatic / Adaptive Narrow - iLock	Automatic / Adaptive Narrow - iLock
Noise Reduction	Noise Management SoundSmoothing Directional	Noise Management SoundSmoothing Directional	Noise Management SoundSmoothing Directional	Noise Management SoundSmoothing Directional	Noise Management SoundSmoothing	Noise Management SoundSmoothing
Wind Noise Reduction	Standard	Standard Binaural	Standard	Standard Binaural	Standard	Standard
Reverb Reducer	•	•	-	-	-	-
Bandwidth: Extension / Compression	• / •	• / •	- / •	- / •	- / •	- / •
Music Enhancer: Live / Recording / Acting	• / • / •	• / • / •	• / - / -	• / - / -	-	-
Tinnitus	Sound Therapy Notch Therapy	Sound Therapy Notch Therapy	Sound Therapy Notch Therapy	Sound Therapy Notch Therapy	Sound Therapy Notch Therapy	Sound Therapy Notch Therapy
Xphone	•	•	•	•	•	•
Acclimatization/Data Logging	Premium / •	Premium / •	Premium / •	Premium / •	High / •	High / •

Accessories & Apps		
Smart Mic	-	•
Smart Transmitter 2,4	-	•
Smart Key	•	•
CROS inoX 8C	•	-
CROS RIC 8C	-	•
Rexton App	•	•



FITTING RANGES

SLIM-RIC HEARING AIDS



STYLELINE





STYLELINE BT



Battery size

Rechargeable battery (Li-Ion)

Rechargeable battery (Li-Ion)

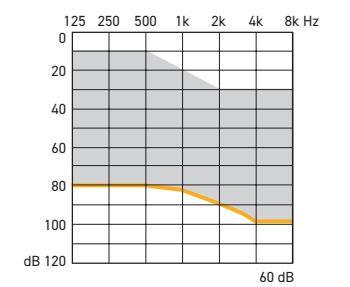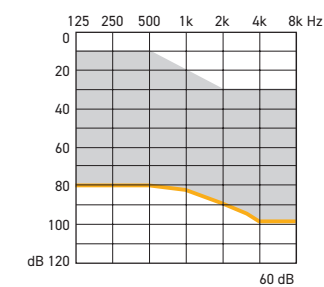

STYLELINE



Performance Level 80 | 60 | 40 | 30 | 20



Highlights


- › Lithium-ion battery
- › Smart Demo device available
- › Programming via: ConnexxAir

INNOVATIVE CHARGING CONCEPT
MEETS NEW SLIM-RIC FORM FACTOR
AND MOBILE CHARGING CASE

Rexton StyleLine - proven to overcome the hearing aid stigma in a new, previously unmatched way. Higher consumer approval than any other hearing aid. the perfect approach for first-time wearers. MyCore technology, innovative design and a new mobile charging system that enables wearers to be on the go for up to four days without needing any power connection, this is the true choice for mobile and active wearers. And with Smart Key or Rexton App, the full functionality of StyleLine can be remotely set and superior hearing quality easily enjoyed.

Receiver	M
FOG/OSPL in dB	60/119
Battery runtime	up to 19 h

Colors

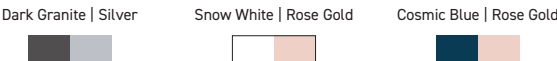

Features / Performance Level	80	60	40	30	20
Channels/Controls/Programs	48/20/6	32/16/6	24/12/6	16/8/4	16/8/4
Soundpro	•	•	•	-	-
My Voice	•	•	•	-	-
Direct Streaming/ Auto Volume	- / -	- / -	- / -	- / -	- / -
Wireless sync	•	•	•	•	•
Directionality	Automatic/Adaptive iOmni Front/Back Left/Right Narrow - Stereo iLock	Automatic/Adaptive iOmni Front/Back Narrow - iLock	Automatic/Adaptive iOmni Narrow - iLock	Automatic/Adaptive iOmni	Automatic/Adaptive iOmni
Noise Reduction	Noise Management SoundSmoothing Directional	Noise Management SoundSmoothing Directional	Noise Management SoundSmoothing	Noise Management SoundSmoothing	Noise Management
Wind Noise Reduction	Standard Binaural	Standard Binaural	Standard	Standard	-
Reverb Reducer	•	-	-	-	-
Bandwidth: Extension / Compression	• / •	- / •	- / •	- / •	- / •
Music Enhancer: Live / Recording / Acting	• / • / •	• / - / -	-	-	-
Tinnitus	Sound Therapy Notch Therapy	Sound Therapy Notch Therapy	Sound Therapy Notch Therapy	Sound Therapy	-
Xphone	•	•	•	-	-
Acclimatization/Data Logging	Premium / •	Premium / •	High / •	High / •	Standard / •
Accessories & Apps					
Mobile charging system	•	•	•	•	•
Smart Key	•	•	•	•	•
Rexton App	•	•	•	•	•

STYLELINE BT



Performance Level 80 | 60 | 40



Highlights


- › Lithium-ion battery
- › Smart Demo device available
- › Programming via: NoahLink Wireless

ENJOY LEADING-EDGE MOBILITY
AND FULL STREAMING WITH
REXTON SLIM-RIC DEVICE

StyleLine BT offers the same reliable feature set as StyleLine, such as reliable MyCore technology and the truly innovative mobile charging system. On top, StyleLine BT integrates direct audio streaming via Bluetooth to iOS devices or other Bluetooth enabled Android systems (via Smart Mic)

Receiver	M
FOG / OSPL in dB	60/119
Battery runtime without streaming	up to 19 h
Battery runtime incl. 5h streaming	up to 16 h

Colors

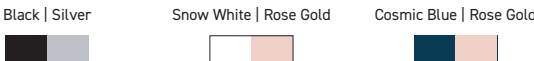

Features / Performance Level	80	60	40
Channels/Controls/Programs	48/20/6	32/16/6	24/12/6
Soundpro	•	•	•
My Voice	•	•	•
Direct Streaming/ Auto Volume	Made for iPhone, Android via TV Transmitter & Smart Mic / Auto Volume	Made for iPhone, Android via TV Transmitter & Smart Mic / Auto Volume	Made for iPhone, Android via TV Transmitter & Smart Mic / Auto Volume
Wireless sync	•	•	•
Directionality	Automatic/Adaptive iOmni Front/Back Left/Right Narrow - Stereo iLock	Automatic/Adaptive iOmni Front/Back Narrow - iLock	Automatic/Adaptive iOmni Narrow - iLock
Noise Reduction	Noise Management SoundSmoothing Directional	Noise Management SoundSmoothing Directional	Noise Management SoundSmoothing
Wind Noise Reduction	Standard Binaural	Standard Binaural	Standard
Reverb Reducer	•	-	-
Bandwidth: Extension / Compression	• / •	- / •	- / •
Music Enhancer: Live / Recording / Acting	• / • / •	• / - / -	-
Tinnitus	Sound Therapy Notch Therapy	Sound Therapy Notch Therapy	Sound Therapy Notch Therapy
Xphone	•	•	•
Acclimatization/Data Logging	Premium / •	Premium / •	High / •
Accessories & Apps			
Mobile charging system	•	•	•
Smart Mic	•	•	•
Smart Transmitter 2,4	•	•	•
Smart Key	•	•	•
Rexton App	•	•	•



FITTING RANGES

MYCORE RIC HEARING AIDS



EMERALD XS



EMERALD S


















STELLAR RIC



EMERALD M



Battery size	10	312	Rechargeable battery (Li-Ion)	13
--------------	----	-----	-------------------------------	----

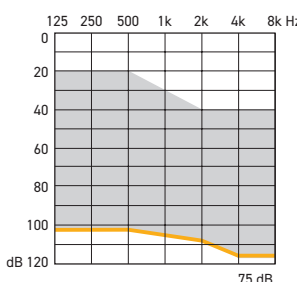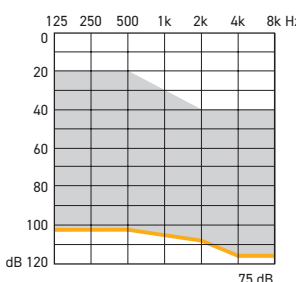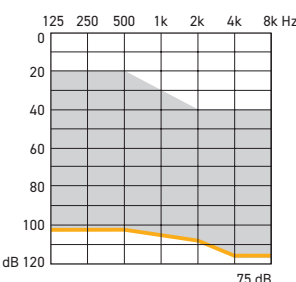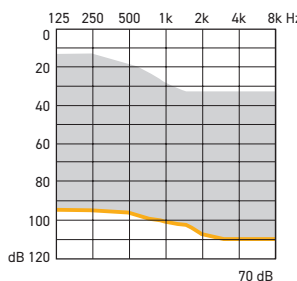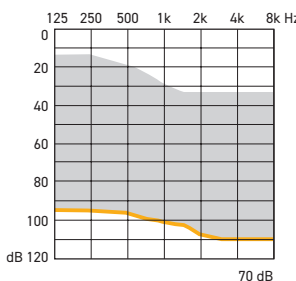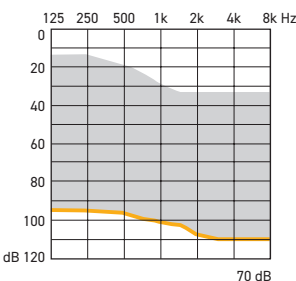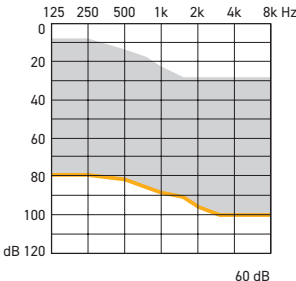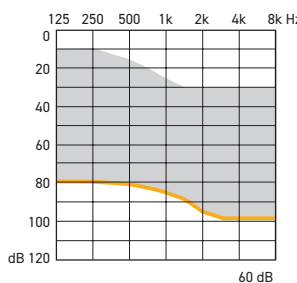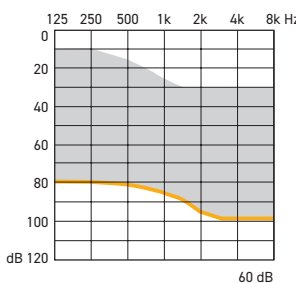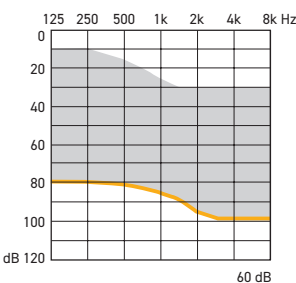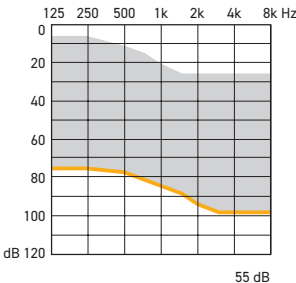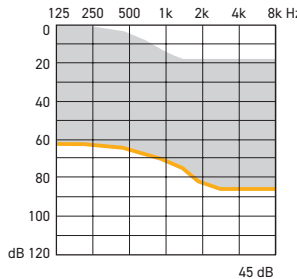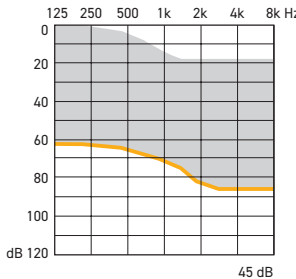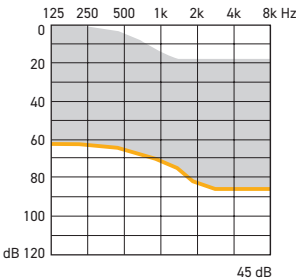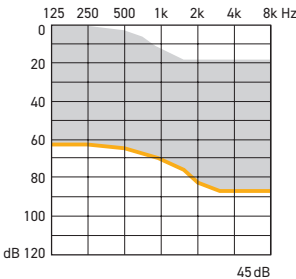

EMERALD XS



Performance Level 80 | 60 | 40 | 30 | 20

OUR SMALLEST RIC WITH NEAR INVISIBILITY BEHIND THE EAR

The Emerald XS is the ideal combination of small size and advanced functionality. Thanks to its durable design, you can enjoy the sounds of nature and music discreetly, no matter your lifestyle or degree or hearing loss. Perfect for first-time users demanding high discretion. Smart Key or Rexton App enable wearers to change volume, programs, and more.



Highlights

- › Battery size 10
- › Mini Receivers 2.0
- › Smart Demo device available
- › Programming via:
- › HiPro

Receiver	S	M	P
FOG/OSPL in dB	45/108	55/113	60/118
Battery current drain	1.2 mA	1.3 mA	1.2 mA

Colors



Features / Performance Level	80	60	40	30	New 20
Channels/Controls/Programs	48/20/6	32/16/6	24/12/6	16/8/4	16/8/4
Soundpro	•	•	•	-	-
My Voice	-	-	-	-	-
Direct Streaming/ Auto Volume	- / -	- / -	- / -	- / -	- / -
Wireless sync	-	-	-	-	-
Directionality	Automatic / Adaptive iOmni Front & Back	Automatic / Adaptive iOmni Front & Back	Automatic / Adaptive iOmni	Automatic / Adaptive iOmni	Automatic / Adaptive iOmni
Noise Reduction	Noise Management SoundSmoothing Directional	Noise Management SoundSmoothing Directional	Noise Management SoundSmoothing	Noise Management SoundSmoothing	Noise Management
Wind Noise Reduction	Standard	Standard	Standard	Standard	-
Reverb Reducer	•	-	-	-	-
Bandwidth: Extension / Compression	• / •	- / •	- / •	- / •	- / •
Music Enhancer: Live / Recording / Acting	• / • / •	• / - / -	-	-	-
Tinnitus	Sound Therapy Notch Therapy	Sound Therapy Notch Therapy	Sound Therapy Notch Therapy	Sound Therapy	-
Xphone	-	-	-	-	-
Acclimatization/Data Logging	Premium / •	Premium / •	High / •	High / •	Standard / •
Accessories & Apps					
Smart Key	•	•	•	•	•
Rexton App	•	•	•	•	•

EMERALD S



Performance Level 80 | 60 | 40 | 30 | 20

A PERFECT COMBINATION OF SMALL SIZE AND PERFORMANCE

This RIC 312 boasts direct audio streaming from iPhone via Bluetooth® and streaming with most Bluetooth® - enabled phones via Smart Mic. Emerald S is small yet fully featured. Its variety of accessories offers full control and usability, e.g. TVs and other Bluetooth® enabled devices via Smart Transmitter 2,4. My Voice ensures a wearer's natural own voice experience by processing it differently than other voices and sounds.



Highlights

- › Battery size 312
- › Mini Receivers 2.0
- › Smart Demo device available
- › Programming via:
- › NoahLink Wireless
- › HiPro

Receiver	S	M	P	HP
FOG/OSPL in dB	45/108	60/119	70/124	75/130
Battery current drain	1.2 mA	1.4 mA	1.3 mA	1.3 mA

Colors



Features / Performance Level	80	60	40	30	New 20
Channels/Controls/Programs	48/20/6	32/16/6	24/12/6	16/8/4	16/8/4
Soundpro	•	•	•	-	-
My Voice	•	•	•	-	-
Direct Streaming/ Auto Volume	Made for iPhone, Android via TV Transmitter & Smart Mic / Auto Volume	Made for iPhone, Android via TV Transmitter & Smart Mic / Auto Volume	Made for iPhone, Android via TV Transmitter & Smart Mic / Auto Volume	Made for iPhone, Android via TV Transmitter & Smart Mic / Auto Volume	Made for iPhone, Android via TV Transmitter & Smart Mic / Auto Volume
Wireless sync	•	•	•	•	•
Directionality	Automatic / Adaptive iOmni Front & Back Left & Right Narrow - Stereo iLock	Automatic / Adaptive iOmni Front & Back Narrow - iLock	Automatic / Adaptive iOmni Narrow - iLock	Automatic / Adaptive iOmni	Automatic / Adaptive iOmni
Noise Reduction	Noise Management SoundSmoothing Directional	Noise Management SoundSmoothing Directional	Noise Management SoundSmoothing	Noise Management SoundSmoothing	Noise Management
Wind Noise Reduction	Standard Binaural	Standard Binaural	Standard	Standard	-
Reverb Reducer	•	-	-	-	-
Bandwidth: Extension / Compression	• / •	- / •	- / •	- / •	- / •
Music Enhancer: Live / Recording / Acting	• / • / •	• / - / -	-	-	-
Tinnitus	Sound Therapy Notch Therapy	Sound Therapy Notch Therapy	Sound Therapy Notch Therapy	Sound Therapy	-
Xphone	•	•	•	-	-
Acclimatization/Data Logging	Premium / •	Premium / •	High / •	High / •	Standard / •
Accessories & Apps					
Smart Mic	•	•	•	•	•
Smart Transmitter 2,4	•	•	•	•	•
Smart Key	•	•	•	•	•
CROS RIC 8C	•	•	•	-	-
Rexton App	•	•	•	•	•

STELLAR RIC



Performance Level 80 | 60 | 40 | 30 | 20



Highlights

- › Lithium-ion battery
- › Mini Receivers 2.0
- › Smart Demo device available
- › Programming via: NoahLink Wireless

LITHIUM-ION RECHARGEABLE RIC THAT PROVIDES MYCORE HEARING THAT LASTS ALL DAY

Stellar RIC is our lithium-ion rechargeable RIC that does not require any battery changes. Stellar RIC is suitable for up to severe hearing loss, and with an IP68 rating it can stand up to moisture, sweat and dirt. Contactless charging for easy, automatic charging that stops when fully charged. My Voice ensures a wearer’s natural own voice experience.

Receiver	S	M	P	HP
FOG/OSPL in dB	45/108	60/119	70/124	75/130
Batt. runtime w/o streaming	up to 21 h	up to 21 h	up to 21 h	up to 21 h
Batt. runtime w. 5h streaming	up to 19 h	up to 19 h	up to 19 h	up to 19 h

Colors

Beige

Deep Brown

Dark Granite

Granite

Gray

Pearl White

Sandy Brown

Silver

Features / Performance Level	80	60	40	30	20
Channels/Controls/Programs	48/20/6	32/16/6	24/12/6	16/8/4	16/8/4
Soundpro	•	•	•	-	-
My Voice	•	•	•	-	-
Direct Streaming/ Auto Volume	Made for iPhone, Android via TV Transmitter & Smart Mic / Auto Volume	Made for iPhone, Android via TV Transmitter & Smart Mic / Auto Volume	Made for iPhone, Android via TV Transmitter & Smart Mic / Auto Volume	Made for iPhone, Android via TV Transmitter & Smart Mic / Auto Volume	Made for iPhone, Android via TV Transmitter & Smart Mic / Auto Volume
Wireless sync	•	•	•	•	•
Directionality	Automatic / Adaptive iOmni Front/Back Left/Right Narrow - Stereo iLock	Automatic / Adaptive iOmni Front/ Back Narrow - iLock	Automatic / Adaptive iOmni Narrow - iLock	Automatic / Adaptive iOmni	Automatic / Adaptive iOmni
Noise Reduction	Noise Management SoundSmoothing Directional	Noise Management SoundSmoothing Directional	Noise Management SoundSmoothing	Noise Management SoundSmoothing	Noise Management
Wind Noise Reduction	Standard Binaural	Standard Binaural	Standard	Standard	-
Reverb Reducer	•	-	-	-	-
Bandwidth: Extension / Compression	• / •	- / •	- / •	- / •	- / •
Music Enhancer: Live / Recording / Acting	• / • / •	• / - / -	-	-	-
Tinnitus	Sound Therapy Notch Therapy	Sound Therapy Notch Therapy	Sound Therapy Notch Therapy	Sound Therapy	-
Xphone	•	•	•	-	-
Acclimatization/Data Logging	Premium / •	Premium / •	High / •	High / •	Standard / •

Accessories & Apps					
Charging Station	•	•	•	•	•
Smart Li-Ion Power	•	•	•	•	•
Smart Mic	•	•	•	•	•
Smart Transmitter 2.4	•	•	•	•	•
Smart Key	•	•	•	•	•
CROS RIC 8C	•	•	•	-	-
CROS Li RIC 8C	•	•	•	-	-
Rexton App	•	•	•	-	-

EMERALD M



Performance Level 80 | 60 | 40 | 30 | 20



Highlights

- › Battery size 13
- › Mini Receivers 2.0
- › Smart Demo device available
- › Programming via: NoahLink Wireless HiPro

GREAT NEW FEATURES AND BENEFITS IN A SMALL HEARING AID

The Emerald M is the ideal combination of small size and long battery lifetime. With up to 75dB of gain, Emerald M fits for even severe hearing loss. The size 13 battery and low power consumption of the hearing aid ensure longest running time. Direct connectivity to iPhones, TVs and other Bluetooth® enabled devices via Smart Transmitter 2.4 and Smart Mic. My Voice ensures a wearer’s natural own voice experience by processing it differently than other voices and sounds. And the optional battery door with integrated T-Coil ensures optimal understanding in locations equipped with loop system.

Receiver	S	M	P	HP
FOG/OSPL in dB	45/108	60/119	70/124	75/130
Battery current drain	1.2 mA	1.4 mA	1.3 mA	1.3 mA

Colors

Beige

Deep Brown

Dark Granite

Granite

Gray

Pearl White

Sandy Brown

Silver

Features / Performance Level	80	60	40	30	20
Channels/Controls/Programs	48/20/6	32/16/6	24/12/6	16/8/4	16/8/4
Soundpro	•	•	•	-	-
My Voice	•	•	•	-	-
Direct Streaming/ Auto Volume	Made for iPhone, Android via TV Transmitter & Smart Mic / Auto Volume	Made for iPhone, Android via TV Transmitter & Smart Mic / Auto Volume	Made for iPhone, Android via TV Transmitter & Smart Mic / Auto Volume	Made for iPhone, Android via TV Transmitter & Smart Mic / Auto Volume	Made for iPhone, Android via TV Transmitter & Smart Mic / Auto Volume
Wireless sync	•	•	•	•	•
Directionality	Automatic / Adaptive iOmni Front & Back Left & Right Narrow - Stereo iLock	Automatic / Adaptive iOmni Front & Back Narrow - iLock	Automatic / Adaptive iOmni Narrow - iLock	Automatic / Adaptive iOmni	Automatic / Adaptive iOmni
Noise Reduction	Noise Management SoundSmoothing Directional	Noise Management SoundSmoothing Directional	Noise Management SoundSmoothing	Noise Management SoundSmoothing	Noise Management
Wind Noise Reduction	Standard Binaural	Standard Binaural	Standard	Standard	-
Reverb Reducer	•	-	-	-	-
Bandwidth: Extension / Compression	• / •	- / •	- / •	- / •	- / •
Music Enhancer: Live / Recording / Acting	• / • / •	• / - / -	-	-	-
Tinnitus	Sound Therapy Notch Therapy	Sound Therapy Notch Therapy	Sound Therapy Notch Therapy	Sound Therapy	-
Xphone	•	•	•	-	-
Acclimatization/Data Logging	Premium / •	Premium / •	High / •	High / •	Standard / •

Accessories & Apps					
Smart Mic	•	•	•	•	•
Smart Transmitter 2,4	•	•	•	•	•
Smart Key	•	•	•	•	•
CROS RIC 8C	•	•	•	-	-
Rexton App	•	•	•	-	-

FITTING RANGES

MYCORE BTE HEARING AIDS









MOSAIC M

STELLAR M

MOSAIC P



Battery size	13	Rechargeable battery (Li-Ion)	13
--------------	----	-------------------------------	----

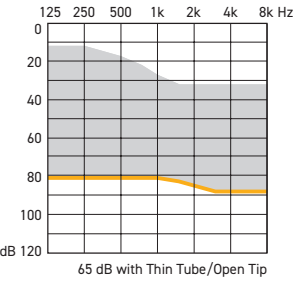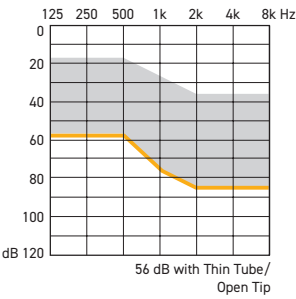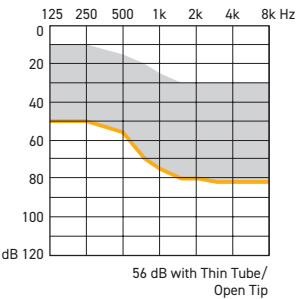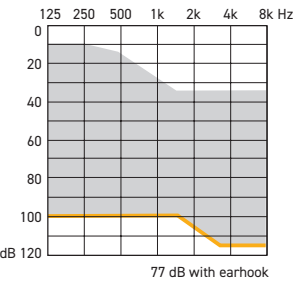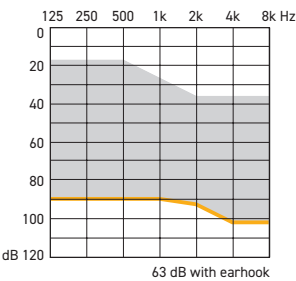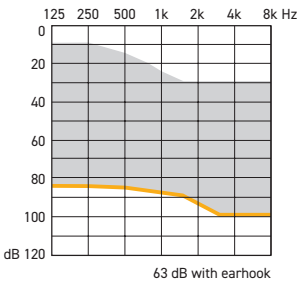

MOSAIC M



Performance Level 80 | 60 | 40 | 30 | 20

THE RELIABLE BTE THAT FITS ANY EAR

The Mosaic M is the ideal combination of small size and long battery lifetime. With up to 63 dB of gain, Mosaic M fits up to severe hearing loss. The size 13 battery and low power consumption ensures longest running time. It offers direct connectivity to iPhones, TVs and other Bluetooth® enabled devices. My Voice ensures a wearer's natural own voice experience by processing it differently than other voices and sounds. Mosaic M is also available with an optional telecoil, integrated in the battery door.



Highlights

- > Battery size 13
- > Smart Demo device available
- > Optional integrated telecoil
- > Programming via:
 - > NoahLink Wireless
 - > HiPro

	Ear hook	ThinTube
FOG/OSPL in dB	63/134	56/126
Battery current drain	1.4 mA	2.0 mA

Colors

Beige

Deep Brown

Dark Granite

Granite

Gray

Pearl White

Sandy Brown

Silver

Features / Performance Level	80	60	40	30 New	20 New
Channels/Controls/Programs	48/20/6	32/16/6	24/12/6	16/8/4	16/8/4
Soundpro	•	•	•	-	-
My Voice	•	•	•	-	-
Direct Streaming/ Auto Volume	Made for iPhone, Android via TV Transmitter & Smart Mic / Auto Volume	Made for iPhone, Android via TV Transmitter & Smart Mic / Auto Volume	Made for iPhone, Android via TV Transmitter & Smart Mic / Auto Volume	Made for iPhone, Android via TV Transmitter & Smart Mic / Auto Volume	Made for iPhone, Android via TV Transmitter & Smart Mic / Auto Volume
Wireless sync	•	•	•	•	•
Directionality	Automatic / Adaptive iOmni Front & Back Left & Right Narrow - Stereo iLock	Automatic / Adaptive iOmni Front & Back Narrow - iLock	Automatic / Adaptive iOmni Narrow - iLock	Automatic / Adaptive iOmni	Automatic / Adaptive iOmni
Noise Reduction	Noise Management SoundSmoothing Directional	Noise Management SoundSmoothing Directional	Noise Management SoundSmoothing	Noise Management SoundSmoothing	Noise Management
Wind Noise Reduction	Standard Binaural	Standard Binaural	Standard	Standard	-
Reverb Reducer	•	-	-	-	-
Bandwidth: Extension / Compression	• / •	- / •	- / •	- / •	- / •
Music Enhancer: Live / Recording / Acting	• / • / •	• / - / -	-	-	-
Tinnitus	Sound Therapy Notch Therapy	Sound Therapy Notch Therapy	Sound Therapy Notch Therapy	Sound Therapy	-
Xphone	•	•	•	-	-
Acclimatization/Data Logging	Premium / •	Premium / •	High / •	High / •	Standard / •
Accessories & Apps					
Smart Mic	•	•	•	•	•
Smart Transmitter 2,4	•	•	•	•	•
Smart Key	•	•	•	•	•
CROS RIC 8C	•	•	•	-	-
Rexton App	•	•	•	•	•

STELLAR M



Performance Level 80 | 60 | 40 | 30 | 20



Highlights

- › Lithium-ion battery
- › Smart Demo device available
- › Programming via: NoahLink Wireless

COMPLETE MYCORE FEATURE SET IN EASY-TO-HANDLE BTE WITH LITHIUM-ION RECHARGEABILITY

With Stellar M, Rexton offers the MyCore feature set including My Voice and direct streaming via Bluetooth® in a robust, easy to handle BTE hearing instrument. Contactless charging for easy, automatic charging that stops when fully charged. Direct connectivity to iPhones, TVs and other Bluetooth®-enabled devices.

	Ear hook	ThinTube
FOG/OSPL in dB	63/134	56/126
Battery runtime without streaming	up to 21 h	up to 21 h
Battery runtime incl. 5h streaming	up to 19 h	up to 19 h

Colors

Beige

Deep Brown

Dark Granite

Granite

Gray

Pearl White

Sandy Brown

Silver

Features / Performance Level	80	60	40	30	20
Channels/Controls/Programs	48/20/6	32/16/6	24/12/6	16/8/4	16/8/4
Soundpro	•	•	•	-	-
My Voice	•	•	•	-	-
Direct Streaming/ Auto Volume	Made for iPhone, Android via TV Transmitter & Smart Mic / Auto Volume	Made for iPhone, Android via TV Transmitter & Smart Mic / Auto Volume	Made for iPhone, Android via TV Transmitter & Smart Mic / Auto Volume	Made for iPhone, Android via TV Transmitter & Smart Mic / Auto Volume	Made for iPhone, Android via TV Transmitter & Smart Mic / Auto Volume
Wireless sync	•	•	•	•	•
Directionality	Automatic / Adaptive iOmni Front & Back Left & Right Narrow - Stereo iLock	Automatic / Adaptive iOmni Front & Back Narrow - iLock	Automatic / Adaptive iOmni Narrow - iLock	Automatic / Adaptive iOmni	Automatic / Adaptive iOmni
Noise Reduction	Noise Management SoundSmoothing Directional	Noise Management SoundSmoothing Directional	Noise Management SoundSmoothing	Noise Management SoundSmoothing	Noise Management
Wind Noise Reduction	Standard Binaural	Standard Binaural	Standard	Standard	-
Reverb Reducer	•	-	-	-	-
Bandwidth: Extension / Compression	• / •	- / •	- / •	- / •	- / •
Music Enhancer: Live / Recording / Acting	• / • / •	• / - / -	-	-	-
Tinnitus	Sound Therapy Notch Therapy	Sound Therapy Notch Therapy	Sound Therapy Notch Therapy	Sound Therapy	-
Xphone	•	•	•	-	-
Acclimatization/Data Logging	Premium / •	Premium / •	High / •	High / •	Standard / •

Accessories & Apps					
Connexx charging station	•	•	•	•	•
Smart Li-Ion Power	•	•	•	•	•
Smart Mic	•	•	•	•	•
Smart Transmitter 2,4	•	•	•	•	•
Smart Key	•	•	•	•	•
CROS RIC 8C	•	•	•	-	-
CROS Li RIC 8C	•	•	•	-	-
Rexton App	•	•	•	•	•

MOSAIC P



Performance Level 80 | 60 | 40 | 30 | 20



Highlights

- › Battery size 13
- › Smart Demo device available
- › Programming via: NoahLink Wireless HiPro

MOSAIC P: POWER - PURE AND SIMPLE

Mosaic P is perfect for patients who need extra amplification. The Mosaic P features a long-lasting battery with low power consumption and the flexibility you need to treat even severe hearing loss. Direct connectivity to iPhones, TVs and other Bluetooth® enabled devices via Smart Transmitter 2.4 and Smart Mic. My Voice ensures a wearer's natural own voice experience by processing it differently than other voices and sounds. Optional battery door with integrated T-Coil ensures optimal understanding in locations equipped with loop system.

	Ear hook	ThinTube
FOG/OSPL in dB	77/135	65/130
Battery current drain	2.0 mA	1.7 mA

Colors

Beige

Deep Brown

Dark Granite

Granite

Gray

Pearl White

Sandy Brown

Silver

Features / Performance Level	80	60	40	30	20
Channels/Controls/Programs	48/20/6	32/16/6	24/12/6	16/8/4	16/8/4
Soundpro	•	•	•	-	-
My Voice	•	•	•	-	-
Direct Streaming/ Auto Volume	Made for iPhone, Android via TV Transmitter & Smart Mic / Auto Volume	Made for iPhone, Android via TV Transmitter & Smart Mic / Auto Volume	Made for iPhone, Android via TV Transmitter & Smart Mic / Auto Volume	Made for iPhone, Android via TV Transmitter & Smart Mic / Auto Volume	Made for iPhone, Android via TV Transmitter & Smart Mic / Auto Volume
Wireless sync	•	•	•	•	•
Directionality	Automatic / Adaptive iOmni Front & Back Left & Right Narrow - Stereo iLock	Automatic / Adaptive iOmni Front & Back Narrow - iLock	Automatic / Adaptive iOmni Narrow - iLock	Automatic / Adaptive iOmni	Automatic / Adaptive iOmni
Noise Reduction	Noise Management SoundSmoothing Directional	Noise Management SoundSmoothing Directional	Noise Management SoundSmoothing	Noise Management SoundSmoothing	Noise Management
Wind Noise Reduction	Standard Binaural	Standard Binaural	Standard	Standard	-
Reverb Reducer	•	-	-	-	-
Bandwidth: Extension / Compression	• / •	- / •	- / •	- / •	- / •
Music Enhancer: Live / Recording / Acting	• / • / •	• / - / -	-	-	-
Tinnitus	Sound Therapy Notch Therapy	Sound Therapy Notch Therapy	Sound Therapy Notch Therapy	Sound Therapy	-
Xphone	•	•	•	-	-
Acclimatization/Data Logging	Premium / •	Premium / •	High / •	High / •	Standard / •

Accessories & Apps					
Smart Mic	•	•	•	•	•
Smart Transmitter 2,4	•	•	•	•	•
Smart Key	•	•	•	•	•
CROS RIC 8C	•	•	•	-	-
Rexton App	•	•	•	•	•



TRUCORE

TRULY CONVINCING TECHNOLOGY



VOICE RANGER

Delivering a more relaxed listening experience



MUSIC ENHANCER

Maximizing listening enjoyment for music lovers



XPHONE

Sending telephone signals automatically to both hearing aids



REVERB REDUCER

Pleasant communication in an echoing room



TRUCORE

CUSTOM HEARING AIDS OVERVIEW



INOX



STERLING CIC













STERLING ITC



STERLING ITE



Battery size	10 / 312	10	312	312
--------------	----------	----	-----	-----

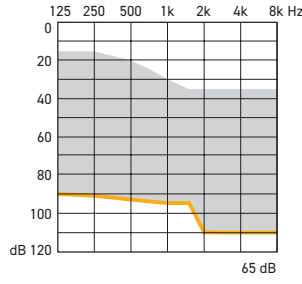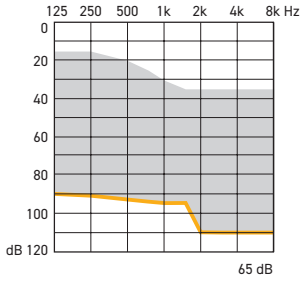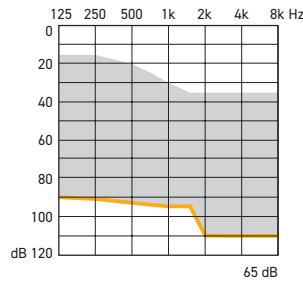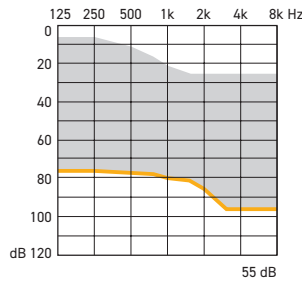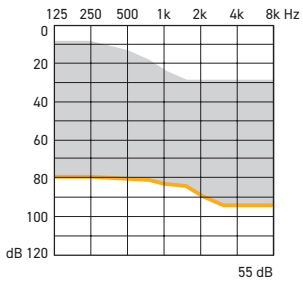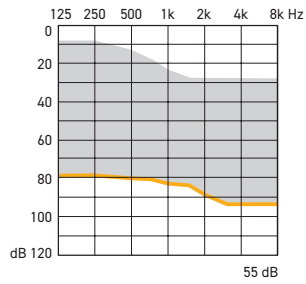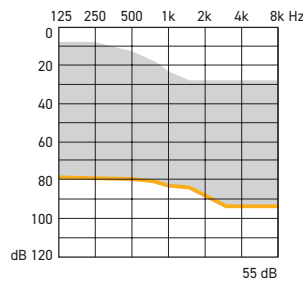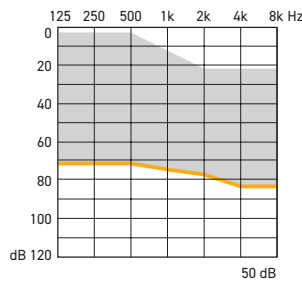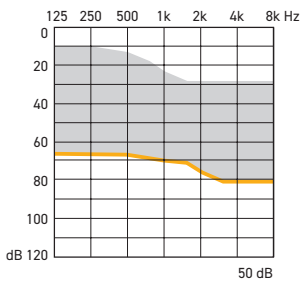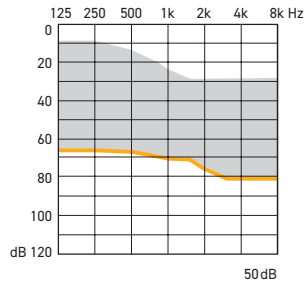



CUSTOM DISCRETION
WITHOUT CUSTOM PRODUCTION

The shape of the inoX housing follows the natural auditory canal geometry and fits for nearly every ear. As discreet as custom CIC's, but without the hassle of custom production. An ideal solution for first-time wearers who favor discretion and convenience above all other. Also available as ITC.

Highlights

- › Battery size 10 (CIC)
- › Battery size 312 (ITC)
- › Programming via:

› Programming cable Flex-Connect

› ConnexxAir (PL80 to 20)

	inoX CIC	inoX ITC
FOG/OSPL in dB	50/113	55/116
Battery current drain	1.1 mA	1.1 mA

Colors

PL 80–20

PL 10

RedBlueMocha

BeigeMocha

Shell

Faceplate

Features / Performance Level	80	60	40	30	20	10
Channels/Controls/Programs	40 / 20 / 6	32 / 16 / 6	24 / 12 / 6	16 / 8 / 4	16 / 8 / 4	12 / 6 / 4
Smart Automatic Equalizer	•	•	•	-	-	-
Wireless Sync	•	•	•	•	•	-
Directionality	Automatic / Adaptive bin. Single-mic Directionality	-	-	-	-	-
Noise Reduction	Noise Management SoundSmoothing Directional	Noise Management SoundSmoothing	Noise Management SoundSmoothing	Noise Management SoundSmoothing	Noise Management	Noise Management
Wind Noise Cancellation	Standard	Standard	Standard	Standard	-	-
Reverb Reducer	•	-	-	-	-	-
Bandwidth: Extension / Compression	• / •	- / •	- / •	- / •	- / •	- / -
Music Enhancer: Live/Recording/Acting	• / • / •	• / - / -	• / - / -	-	-	-
Tinnitus Function	Sound Therapy	Sound Therapy	Sound Therapy	Sound Therapy	-	-
Xphone	-	-	-	-	-	-
Acclimatization / Data Logging	Premium / •	Premium / •	High / •	High / •	Standard / •	- / •

Accessories & Apps						
Smart Connect	•	•	•	•	•	-
Smart Transmitter	•	•	•	•	•	-
Smart Remote	•	•	•	•	•	-
Smart Key	•	•	•	•	•	•
CROS RIC	•	•	•	-	-	-
Smart Remote App	•	•	•	•	•	•
Smart Connect App	•	•	•	•	•	-
Rexton App	•	•	•	•	•	•



PROVEN DESIGN

This hearing system provides a pleasant sound experience with excellent listening comfort and high amplification power. Wireless functionality up to 65 dB possible. Binaural signal processing ensures easier listening in loud environments. Smart Key or Rexton App enables wearers to discreetly change volume, programs, and more. Sterling ITC is available with an optional telecoil.

Highlights

- › Battery size 10
- › Programming via:

› Programming cable Flex-Connect

› ConnexxAir (80 to 20 wireless)

	50 dB	55 dB	65 dB
FOG/OSPL in dB	50/113	55/118	65/124
Battery current drain	0.8 mA	0.8 mA	1.0 mA

Colors

BeigeBrownTanMochaDark Brown

Shell

Faceplate

Features / Performance Level	80	60	40	30	20	10
Channels/Controls/Programs	40 / 20 / 6	32 / 16 / 6	24 / 12 / 6	16 / 8 / 4	16 / 8 / 4	12 / 6 / 4
Smart Automatic Equalizer	•	•	•	-	-	-
Wireless Sync	•	•	•	•	•	-
Directionality	Automatic / Adaptive bin. Single-mic Directionality	-	-	-	-	-
Noise Reduction	Noise Management SoundSmoothing Directional	Noise Management SoundSmoothing	Noise Management SoundSmoothing	Noise Management SoundSmoothing	Noise Management	Noise Management
Wind Noise Cancellation	Standard	Standard	Standard	Standard	-	-
Reverb Reducer	•	-	-	-	-	-
Bandwidth: Extension / Compression	• / •	- / •	- / •	- / •	- / •	- / -
Music Enhancer: Live/Recording/Acting	• / • / •	• / - / -	• / - / -	-	-	-
Tinnitus Function	Sound Therapy	Sound Therapy	Sound Therapy	Sound Therapy	-	-
Xphone	-	-	-	-	-	-
Acclimatization / Data Logging	Premium / •	Premium / •	High / •	High / •	Standard / •	- / •

Accessories & Apps						
Smart Connect	•	•	•	•	•	-
Smart Transmitter	•	•	•	•	•	-
Smart Remote	•	•	•	•	•	-
Smart Key	•	•	•	•	•	•
CROS RIC	•	•	•	-	-	-
Smart Remote App	•	•	•	•	•	•
Smart Connect App	•	•	•	•	•	-
Rexton App	•	•	•	•	•	•

STERLING ITC



Performance Level 80 | 60 | 40 | 30 | 20 | 10



MODERN AND COMFORTABLE

The high-performance model with programmable push button. Wireless functionality up to 65 dB possible. Binaural signal processing ensures easier listening in loud environments. Smart Key or the Rexton App enable wearers to discreetly volume, programs and more.

Highlights

- › Battery size 312 and 13
- › Optional telecoil
- › Programming via:

› Programming cable Flex-Connect

› ConnexxAir (80 to 20 wireless)

	50 dB	55 dB	65 dB
FOG/OSPL in dB	50/113	55/118	65/124
Battery current drain	0.8 mA	0.9 mA	1.0 mA

Colors

	Beige	Brown	Tan	Mocha	Dark Brown
Shell					
Faceplate					

Features / Performance Level	80	60	40	30	20	10
Channels/Controls/Programs	40 / 20 / 6	32 / 16 / 6	24 / 12 / 6	16 / 8 / 4	16 / 8 / 4	12 / 6 / 4
Smart Automatic Equalizer	•	•	•	-	-	-
Wireless Sync	•	•	•	•	•	-
Directionality	•	•	•	•	•	-
Noise Reduction	Noise Management SoundSmoothing Directional	Noise Management SoundSmoothing Directional	Noise Management SoundSmoothing	Noise Management SoundSmoothing	Noise Management	Noise Management
Wind Noise Cancellation	Premium	High	Standard	Standard	-	-
Reverb Reducer	•	-	-	-	-	-
Bandwidth: Extension / Compression	• / •	- / •	- / •	- / •	- / •	- / -
Music Enhancer: Live/Recording/Acting	• / • / •	• / - / -	• / - / -	-	-	-
Tinnitus Function	Sound Therapy	Sound Therapy	Sound Therapy	Sound Therapy	-	-
Xphone	•	•	•	-	-	-
Acclimatization / Data Logging	Premium / •	Premium / •	High / •	High / •	Standard / •	- / •
Accessories & Apps						
Smart Connect	•	•	•	•	•	-
Smart Transmitter	•	•	•	•	•	-
Smart Remote	•	•	•	•	•	-
Smart Key	•	•	•	•	•	•
CROS RIC	•	•	•	-	-	-
Smart Remote App	•	•	•	•	•	•
Smart Connect App	•	•	•	•	•	-
Rexton App	•	•	•	•	•	•

STERLING ITE



Performance Level 80 | 60 | 40 | 30 | 20 | 10



INNOVATIVE AND UNIQUE

Our first custom model with 2 microphones embedded in the faceplate delivers a new dimension of sound quality. And connected with the Rexton app, Sterling ITE is a most discreet hearing aid and more. Sterling ITE is available with an optional telecoil.

Highlights

- › Battery size 312 and 13
- › Optional telecoil
- › Programming via:

› Programming cable Flex-Connect

› ConnexxAir (80 to 20 wireless)

	55 dB	65 dB
FOG/OSPL in dB	55/118	65/124
Battery current drain	0.9 mA	1.0 mA

Colors

	Beige	Brown	Tan	Mocha	Dark Brown
Shell					
Faceplate					

Features / Performance Level	80	60	40	30	20	10
Channels/Controls/Programs	40 / 20 / 6	32 / 16 / 6	24 / 12 / 6	16 / 8 / 4	16 / 8 / 4	12 / 6 / 4
Smart Automatic Equalizer	•	•	•	-	-	-
Wireless Sync	•	•	•	•	•	-
Directionality	•	•	•	•	•	-
Noise Reduction	Noise Management SoundSmoothing Directional	Noise Management SoundSmoothing Directional	Noise Management SoundSmoothing	Noise Management SoundSmoothing	Noise Management	Noise Management
Wind Noise Cancellation	Premium	High	Standard	Standard	-	-
Reverb Reducer	•	-	-	-	-	-
Bandwidth: Extension / Compression	• / •	- / •	- / •	- / •	- / •	- / -
Music Enhancer: Live/Recording/Acting	• / • / •	• / - / -	• / - / -	-	-	-
Tinnitus Function	Sound Therapy	Sound Therapy	Sound Therapy	Sound Therapy	-	-
Xphone	•	•	•	-	-	-
Acclimatization / Data Logging	Premium / •	Premium / •	High / •	High / •	Standard / •	- / •
Accessories & Apps						
Smart Connect	•	•	•	•	•	-
Smart Transmitter	•	•	•	•	•	-
Smart Remote	•	•	•	•	•	-
Smart Key	•	•	•	•	•	•
CROS RIC	•	•	•	-	-	-
Smart Remote App	•	•	•	•	•	•
Smart Connect App	•	•	•	•	•	-
Rexton App	•	•	•	•	•	•

TRUCORE RIC HEARING AIDS

EMERALD XS



EMERALD S



EMERALD M



STELLAR LI


















Battery size

10

312 / rechargeable battery

13 / rechargeable battery

Rechargeable battery (Li-Ion)

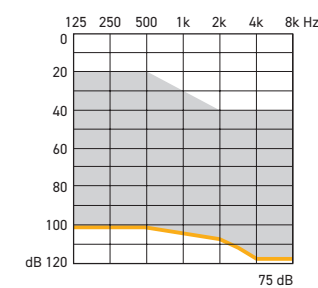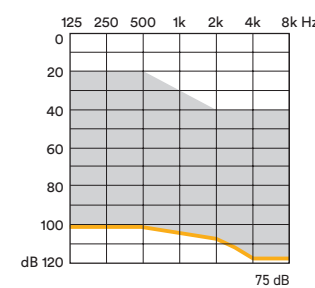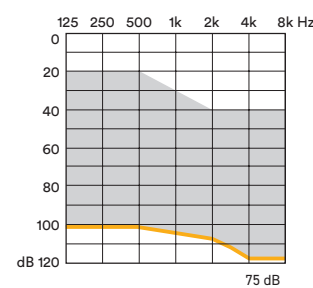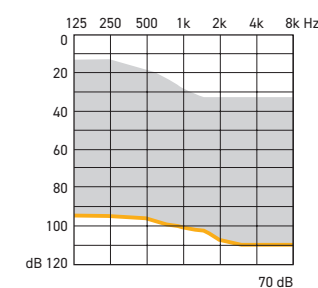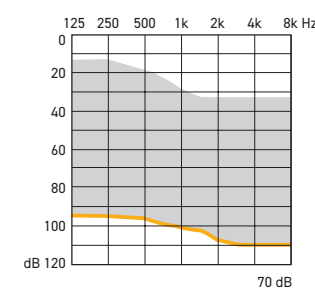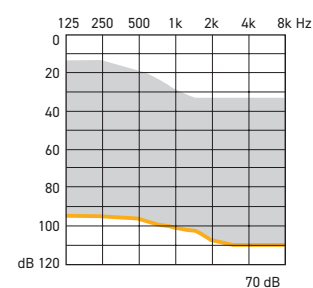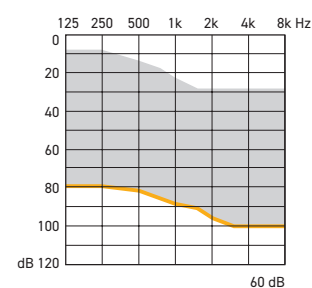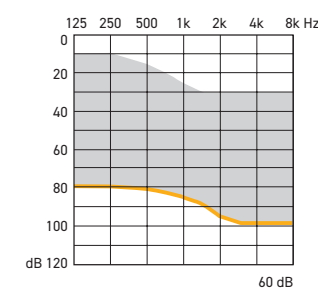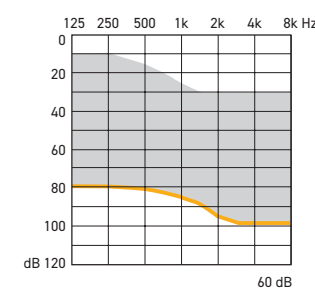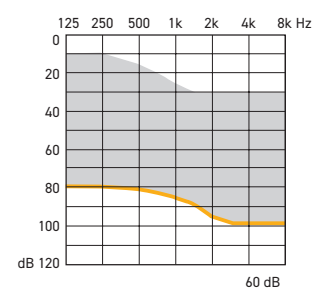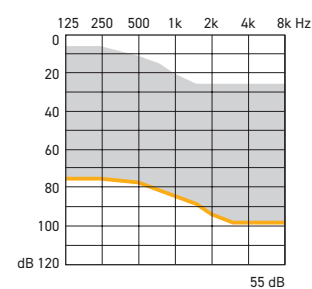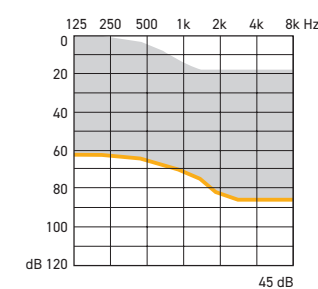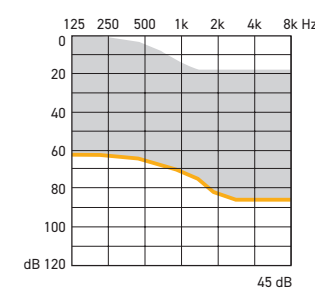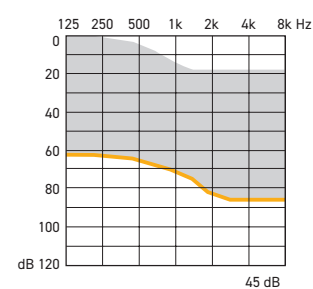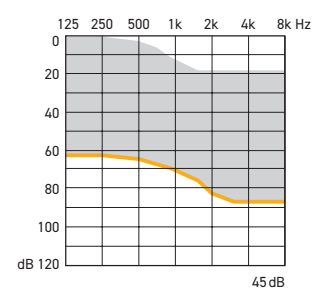

EMERALD XS



Performance Level 80 | 60 | 40 | 30 | 20



SMALL AND UNOBTUSIVE

The Emerald XS hearing aids are equipped with many automatic functions which ensure outstanding hearing comfort. They adapt flexibly to each hearing situation. Perfect for first-time users demanding high discretion. Smart Key or Rexton App enables wearers to discreetly volume, programs and more.

Highlights

- › Battery size 10
- › Mini Receivers 2.0
- › Programming via:
› Programming adaptor 10

Receiver	S	M	P
FOG/OSPL in dB	45/108	55/113	60/118
Battery current drain	0.9 mA	1.0 mA	0.9 mA

Colors
PL 80-40

Beige

Brown

Dark Granite

Granite

Gray

Pearl White

Sandy Brown

Silver

Colors
PL 30-20

Beige

Brown

Dark Granite

Gray

Silver

Features / Performance Level	80	60	40	30	20
Channels/Controls/Programs	40 / 20 / 6	32 / 16 / 6	24 / 12 / 6	16 / 8 / 4	16 / 8 / 4
Smart Automatic Equalizer	•	•	•	-	-
Wireless Sync	-	-	-	-	-
Directionality	Automatic / Adaptive iOmni Front & Back	Automatic / Adaptive iOmni Front & Back	Automatic / Adaptive iOmni	Automatic / Adaptive iOmni	Automatic iOmni
Noise Reduction	Noise Management SoundSmoothing Directional	Noise Management SoundSmoothing Directional	Noise Management SoundSmoothing	Noise Management SoundSmoothing	Noise Management
Wind Noise Cancellation	Standard	Standard	Standard	Standard	-
Reverb Reducer	•	-	-	-	-
Bandwidth: Extension / Compression	• / •	- / •	- / •	- / •	- / •
Music Enhancer: Live/Recording/Acting	• / • / •	• / - / -	• / - / -	-	-
Tinnitus Function	Sound Therapy	Sound Therapy	Sound Therapy	Sound Therapy	-
Xphone	-	-	-	-	-
Acclimatization / Data Logging	Premium / •	Premium / •	High / •	High / •	Standard / •
Accessories & Apps					
Smart Connect	-	-	-	-	-
Smart Transmitter	-	-	-	-	-
Smart Remote	-	-	-	-	-
Smart Key	•	•	•	•	•
CROS RIC	-	-	-	-	-
Smart Remote App	•	•	•	•	•
Rexton App	•	•	•	•	•

EMERALD S



Performance Level 80 | 60 | 40 | 30 | 20 | 10



THE MULTITALENTED DEVICE FOR HIGHEST EXPECTATIONS

The versatile Emerald S is a small, inconspicuous hearing system. Various types of receiver mean that it can be used for a broad range of hearing loss. The Emerald S hearing aids in Performance Levels 80 to 40 can be used with either a standard battery OR rechargeable battery.

Highlights

- › Battery size 312
- › Mini Receivers 2.0
- › Rechargeable battery (80 / 60 /40)
- › Programming via:
› Programming adaptor 312
› ConnexxAir (PL 80-20)

Receiver	S	M	P	HP
FOG/OSPL in dB	45/108	60/119	70/124	75/130
Battery current drain	0.9 mA	1.0 mA	1.0 mA	1.1 mA

Colors
PL 80-40

Beige

Brown

Dark Granite

Granite

Gray

Pearl White

Sandy Brown

Silver

Colors
PL 30-20

Beige

Brown

Dark Granite

Gray

Silver

Features / Performance Level	80	60	40	30	20	10
Channels/Controls/Programs	40 / 20 / 6	32 / 16 / 6	24 / 12 / 6	16 / 8 / 4	16 / 8 / 4	12 / 6 / 4
Smart Automatic Equalizer	•	•	•	-	-	-
Wireless Sync	•	•	•	•	•	-
Directionality	Automatic / Adaptive iOmni Front & Back Left & Right Narrow - iLock	Automatic / Adaptive iOmni Front & Back Narrow - iLock	Automatic / Adaptive iOmni Narrow - iLock	Automatic / Adaptive iOmni	Automatic iOmni	Automatic iOmni
Noise Reduction	Noise Management SoundSmoothing Directional	Noise Management SoundSmoothing Directional	Noise Management SoundSmoothing	Noise Management SoundSmoothing	Noise Management	Noise Management
Wind Noise Cancellation	Standard Binaural	Standard Binaural	Standard	Standard	-	-
Reverb Reducer	•	-	-	-	-	-
Bandwidth: Extension / Compression	• / •	- / •	- / •	- / •	- / •	- / -
Music Enhancer: Live/Recording/Acting	• / • / •	• / - / -	• / - / -	-	-	-
Tinnitus Function	Sound Therapy	Sound Therapy	Sound Therapy	Sound Therapy	-	-
Xphone	•	•	•	-	-	-
Acclimatization / Data Logging	Premium / •	Premium / •	High / •	High / •	Standard / •	- / •
Accessories & Apps						
Smart Connect	•	•	•	•	•	-
Smart Transmitter	•	•	•	•	•	-
Smart Remote	•	•	•	•	•	-
Smart Key	•	•	•	•	•	•
CROS RIC	•	•	•	-	-	-
Smart Remote App	•	•	•	•	•	•
Smart Connect App	•	•	•	•	•	-
Rexton App	•	•	•	•	•	•
Smart Power Charger	•	•	•	-	-	-

EMERALD M



Performance Level 80 | 60 | 40



COMPACT AND INTELLIGENT
PREMIUM SOLUTIONS

The high-performance Emerald M is characterized by its versatility. The housing contains a highly effective hearing system which can alternatively be used with a standard or rechargeable battery. It offers smart connection options. Telephone coil ensures optimal understanding in locations equipped with loop system. The Emerald M hearing aids in Performance Levels 80 to 40 can be used with either a standard battery OR rechargeable battery.

Highlights

- › Battery size 13
- › Mini Receivers 2.0
- › Rechargeable battery
- › Programming via:
- › Programming adaptor 13
- › ConnexxAir

Receiver	S	M	P	HP
FOG/OSPL in dB	45/108	60/119	70/124	75/135
Battery current drain	0.9 mA	1.1 mA	1.0 mA	1.1 mA

Colors



Features / Performance Level	80	60	40
Channels/Controls/Programs	40 / 20 / 6	32 / 16 / 6	24 / 12 / 6
Smart Automatic Equalizer	•	•	•
Wireless Sync	•	•	•
Directionality	Automatic / Adaptive iOmni Front & Back Left & Right Narrow - iLock	Automatic / Adaptive iOmni Front & Back Narrow - iLock	Automatic / Adaptive iOmni Narrow
Noise Reduction	Noise Management SoundSmoothing Directional	Noise Management SoundSmoothing Directional	Noise Management SoundSmoothing Narrow - iLock
Wind Noise Cancellation	Standard Binaural	Standard Binaural	Standard
Reverb Reducer	•	-	-
Bandwidth: Extension / Compression	• / •	- / •	- / •
Music Enhancer: Live/Recording/Acting	• / • / •	• / - / -	• / - / -
Tinnitus Function	Sound Therapy	Sound Therapy	Sound Therapy
Xphone	•	•	•
Acclimatization / Data Logging	Premium / •	Premium / •	High / •

Accessories & Apps			
Smart Connect	•	•	•
Smart Transmitter	•	•	•
Smart Remote	•	•	•
Smart Key	•	•	•
CROS RIC	•	•	•
Smart Remote App	•	•	•
Smart Connect App	•	•	•
Rexton App	•	•	•
Smart Power Charger	•	•	•

STELLAR LI



Performance Level 80 | 60 | 40 | 30



POWER AND ENDURANCE
FOR 24 HOURS

The Stellar Li has an integrated lithium-ion battery which is inductively charged in the charging station. The battery offers an operating time of 24 hours with full connectivity, with only 4 hours charging time. 30 minute fast-charging allows for 7 hours run time.

Highlights


- › Lithium-ion battery
- › Mini Receivers 2.0
- › Programming via:
- › HiPro CS 44
- › ConnexxAir

Receiver	S	M	P	HP
FOG/OSPL in dB	45/108	60/119	70/124	75/130
Battery current drain	0.9 mA	1.0 mA	1.0 mA	1.1 mA

Colors



Colors

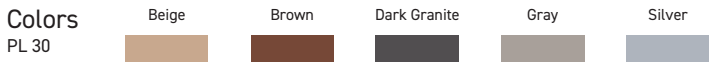

Features / Performance Level	80	60	40	30
Channels/Controls/Programs	40 / 20 / 6	32 / 16 / 6	24 / 12 / 6	16 / 8 / 4
Smart Automatic Equalizer	•	•	•	-
Wireless Sync	•	•	•	•
Directionality	Automatic / Adaptive iOmni Front & Back Left & Right Narrow - iLock	Automatic / Adaptive iOmni Front & Back Narrow - iLock	Automatic / Adaptive iOmni Narrow Narrow - iLock	Automatic / Adaptive iOmni
Noise Reduction	Noise Management SoundSmoothing Directional	Noise Management SoundSmoothing Directional	Noise Management SoundSmoothing	Noise Management SoundSmoothing
Wind Noise Cancellation	Standard Binaural	Standard Binaural	Standard	Standard
Reverb Reducer	•	-	-	-
Bandwidth: Extension / Compression	• / •	- / •	- / •	- / •
Music Enhancer: Live/Recording/Acting	• / • / •	• / - / -	• / - / -	-
Tinnitus Function	Sound Therapy	Sound Therapy	Sound Therapy	Sound Therapy
Xphone	•	•	•	-
Acclimatization / Data Logging	Premium / •	Premium / •	High / •	High / •

Accessories & Apps				
Smart Connect	•	•	•	•
Smart Transmitter	•	•	•	•
Smart Remote	•	•	•	•
Smart Key	•	•	•	•
CROS RIC	•	•	•	-
Smart Remote App	•	•	•	•
Smart Connect App	•	•	•	•
Rexton App	•	•	•	•
Li-ion Charger	•	•	•	•

TRUCORE
TRUCORE BTE HEARING AIDS

MOSAIC S



Battery size 312

MOSAIC M



13 / rechargeable battery

MOSAIC P














13

MOSAIC HP



675

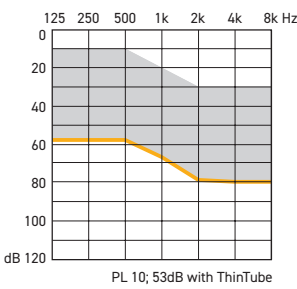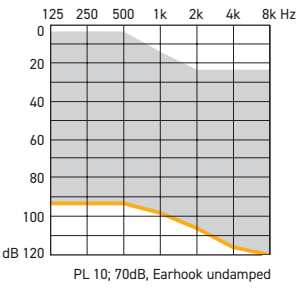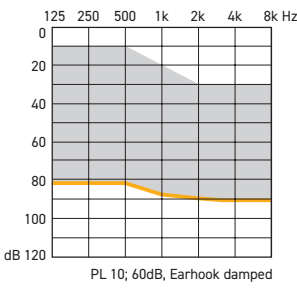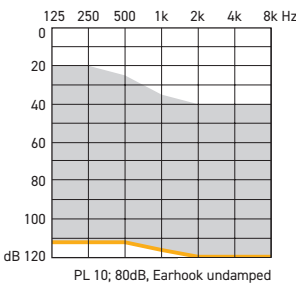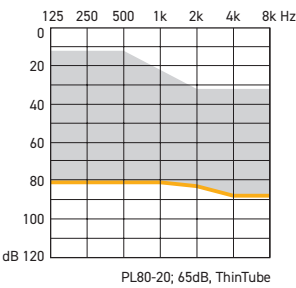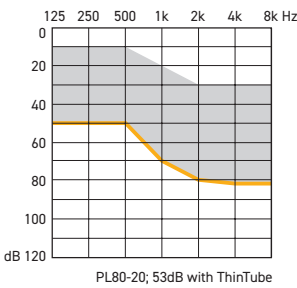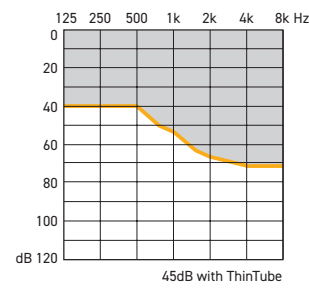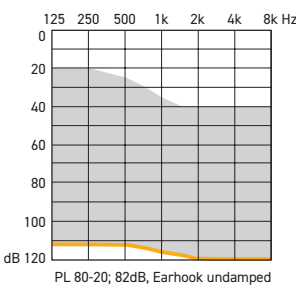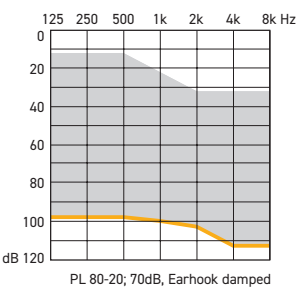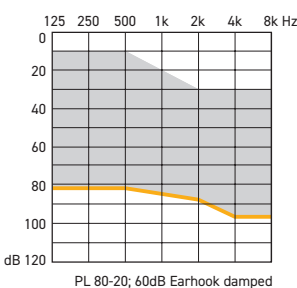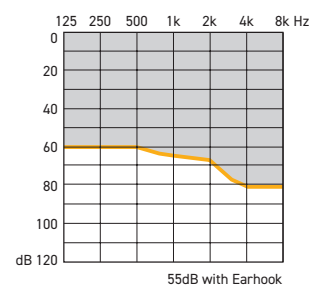

MOSAIC S



Performance Level 40 | 30 | 20 | 10



Highlights

- › Battery size 312
- › Programming via:
 - › ConnexxAir (PL 40-20)
 - › Programming adaptor 312

GRACEFUL AND VARIABLE

Mosaic S - the behind-the-ear hearing system for anyone who particularly values discretion. The Mosaic S is ideally suited for an open-fitting system. This hearing system is a perfect solution for mild to moderate hearing loss.

	Ear hook	ThinTube
F0G/OSPL in dB	55/124	45/124
Battery current drain	0.9 mA	0.9 mA

Colors



Features / Performance Level	40	30	20	10
Channels/Controls/Programs	24 / 12 / 6	16 / 8 / 4	16 / 8 / 4	12 / 6 / 4
Smart Automatic Equalizer	•	-	-	-
Wireless Sync	•	•	•	-
Directionality	Automatic / Adaptive iOmni Narrow - iLock	Automatic / Adaptive iOmni	Automatic iOmni	Automatic iOmni
Noise Reduction	Noise Management SoundSmoothing	Noise Management SoundSmoothing	Noise Management	Noise Management
Wind Noise Cancellation	Standard	Standard	-	-
Reverb Reducer	-	-	-	-
Bandwidth: Extension / Compression	- / •	- / •	- / •	- / -
Music Enhancer: Live/Recording/Acting	• / - / -	-	-	-
Tinnitus Function	Sound Therapy	Sound Therapy	-	-
Xphone	-	-	-	-
Acclimatization / Data Logging	High / •	High / •	Standard / •	- / •

Accessories & Apps				
Smart Connect	•	•	•	-
Smart Transmitter	•	•	•	-
Smart Remote	•	•	•	-
Smart Key	•	•	•	•
Smart Remote App	•	•	•	•
Smart Connect App	•	•	•	-
Rexton App	•	•	•	•

MOSAIC M



Performance Level 80 | 60 | 40 | 30 | 20 | 10



Highlights

- › Battery size 13
- › Optional battery door with Audio Input
- › Rechargeable battery (80 / 60 / 40)
- › Programming via: ConnexxAir (PL 80-20)
- › Programming adaptor 13

UNIVERSAL APPLICATION

The Mosaic M hearing aids are suitable for various degrees of hearing loss and with their many features, they offer wide scope for meeting your customers’ individual needs. Performance Level 80–40 can be used with standard or rechargeable batteries. For open and closed fitting.

	Earhook damped	ThinTube
FOG/OSPL in dB	60/130	53/126 (PL10; 53/125)
Battery current drain	1.0 mA	1.1 mA (PL 10; 1,0 mA)

Colors PL 80–40

Beige

Brown

Dark Granite

Granite

Gray

Pearl White

Sandy Brown

Silver

Colors PL 30-10

Beige

Brown

Dark Granite

Gray

Silver

Features / Performance Level	80	60	40	30	20	10
Channels/Controls/Programs	40 / 20 / 6	32 / 16 / 6	24 / 12 / 6	16 / 8 / 4	16 / 8 / 4	12 / 6 / 4
Smart Automatic Equalizer	•	•	•	-	-	-
Wireless Sync	•	•	•	•	•	-
Directionality	Automatic / Adaptive iOmni Front & Back Left & Right Narrow - iLock	Automatic / Adaptive iOmni Front & Back Narrow - iLock	Automatic / Adaptive iOmni Narrow - iLock	Automatic / Adaptive iOmni	Automatic iOmni	Automatic iOmni
Noise Reduction	Noise Management SoundSmoothing Directional	Noise Management SoundSmoothing Directional	Noise Management SoundSmoothing	Noise Management SoundSmoothing	Noise Management	Noise Management
Wind Noise Cancellation	Standard Binaural	Standard Binaural	Standard	Standard	-	-
Reverb Reducer	•	-	-	-	-	-
Bandwidth: Extension / Compression	• / •	- / •	- / •	- / •	- / •	- / -
Music Enhancer: Live/Recording/Acting	• / • / •	• / - / -	• / - / -	-	-	-
Tinnitus Function	Sound Therapy	Sound Therapy	Sound Therapy	Sound Therapy	-	-
Xphone	•	•	•	-	-	-
Acclimatization / Data Logging	Premium / •	Premium / •	High / •	High / •	Standard / •	- / •

Accessories & Apps						
Smart Connect	•	•	•	•	•	-
Smart Transmitter	•	•	•	•	•	-
Smart Remote	•	•	•	•	•	-
Smart Key	•	•	•	•	•	•
CROS RIC	•	•	•	-	-	-
Smart Remote App	•	•	•	•	•	•
Smart Connect App	•	•	•	•	•	-
Rexton App	•	•	•	•	•	•
Smart Power Charger	•	•	•	-	-	-

MOSAIC P



Performance Level 80 | 60 | 40 | 30 | 20 | 10



Highlights

- › Battery size 13
- › Optional battery door with Audio Input
- › Programming via: ConnexxAir (PL 80-20)
- › Programming adaptor 13 (PL 80-20) and programming shoe (PL 10)

ADVANCED AND CONVINCING

Rexton's Mosaic P models are hearing aids for moderate to severe and severe hearing losses. The Mosaic P model with its rocker switch is very easy to use.

PL 80-20	Earhook damped	Earhook undamped	ThinTube
FOG/OSPL in dB	70/131	77/135	65/129
Battery current drain	1.0 mA	1.6 mA	1.2mA

PL 10	Earhook undamped
FOG/OSPL in dB	70/134
Battery current drain	1.4 mA

Colors PL 80–40

Beige

Brown

Dark Granite

Granite

Gray

Pearl White

Sandy Brown

Silver

Colors PL 30-20

Beige

Deep Brown

Dark Granite

Gray

Silver

Colors PL 10

Beige

Deep Brown

Gray

Features / Performance Level	80	60	40	30	20	10
Channels/Controls/Programs	40 / 20 / 6	32 / 16 / 6	24 / 12 / 6	16 / 8 / 4	16 / 8 / 4	12 / 6 / 4
Smart Automatic Equalizer	•	•	•	-	-	-
Wireless Sync	•	•	•	•	•	-
Directionality	Automatic / Adaptive iOmni Front & Back Left & Right Narrow - iLock	Automatic / Adaptive iOmni Front & Back Narrow - iLock	Automatic / Adaptive iOmni Narrow - iLock	Automatic / Adaptive iOmni	Automatic iOmni	Automatic iOmni
Noise Reduction	Noise Management SoundSmoothing Directional	Noise Management SoundSmoothing Directional	Noise Management SoundSmoothing	Noise Management SoundSmoothing	Noise Management	Noise Management
Wind Noise Cancellation	Standard Binaural	Standard Binaural	Standard	Standard	-	-
Reverb Reducer	•	-	-	-	-	-
Bandwidth: Extension / Compression	• / •	- / •	- / •	- / •	- / •	- / -
Music Enhancer: Live/Recording/Acting	• / • / •	• / - / -	• / - / -	-	-	-
Tinnitus Function	Sound Therapy	Sound Therapy	Sound Therapy	Sound Therapy	-	-
Xphone	•	•	•	-	-	-
Acclimatization / Data Logging	Premium / •	Premium / •	High / •	High / •	Standard / •	- / •

Accessories & Apps						
Smart Connect	•	•	•	•	•	-
Smart Transmitter	•	•	•	•	•	-
Smart Remote	•	•	•	•	•	-
Smart Key	•	•	•	•	•	•
CROS RIC	•	•	•	-	-	-
Smart Remote App	•	•	•	•	•	•
Smart Connect App	•	•	•	•	•	-
Rexton App	•	•	•	•	•	•

MOSAIC HP



Performance Level 80 | 60 | 40 | 30 | 20 | 10



Highlights

- › Battery size 675
- › Optional battery door with Audio Input
- › Child lock
- › Small ear hook
- › Programming via:
 - › ConnexxAir (PL 80-20)
 - › HiPro CS 44 (PL 80-20) and programming shoe (PL 10)

HIGH POWER FOR
ADVANCEDHEARING LOSS

Designed for users who need extra amplification. Mosaic HP offers a long lasting battery with low power consumption and the flexibility you need to treat up to profound hearing losses. Telephone-coil, audio streaming, and useful accessories provide the ultimate flexibility for the wearer, the Rexton App allows for discreet control.

Earhook undamped	
FOG/OSPL in dB	82/140 (PL 10: 80/140)
Battery current drain	2.1 mA (PL10; 2.4 mA)

Colors
PL 80-40

Beige

Brown

Dark Granite

Granite

Gray

Pearl White

Sandy Brown

Silver

Colors
PL 30-20

Beige

Deep Brown

Dark Granite

Gray

Silver

Colors
PL 10

Beige

Deep Brown

Gray

Features / Performance Level	80	60	40	30	20	10
Channels/Controls/Programs	24 / 16 / 6	24 / 16 / 6	24 / 12 / 6	16 / 8 / 4	16 / 8 / 4	12 / 6 / 4
Smart Automatic Equalizer	•	•	•	-	-	-
Wireless Sync	•	•	•	•	•	-
Directionality	Automatic / Adaptive iOmni Front & Back Left & Right Narrow - iLock	Automatic / Adaptive iOmni Front & Back Narrow - iLock	Automatic / Adaptive iOmni Narrow - iLock	Automatic / Adaptive iOmni	Automatic iOmni	Automatic iOmni
Noise Reduction	Noise Management SoundSmoothing Directional	Noise Management SoundSmoothing Directional	Noise Management SoundSmoothing	Noise Management SoundSmoothing	Noise Management	Noise Management
Wind Noise Cancellation	Standard Binaural	Standard Binaural	Standard	Standard	-	-
Reverb Reducer	•	-	-	-	-	-
Bandwidth: Extension / Compression	- / •	- / •	- / •	- / •	- / •	- / -
Music Enhancer: Live/Recording/Acting	• / • / •	• / - / -	• / - / -	-	-	-
Tinnitus Function	Sound Therapy	Sound Therapy	Sound Therapy	Sound Therapy	-	-
Xphone	•	•	•	-	-	-
Acclimatization / Data Logging	Premium / •	Premium / •	High / •	High / •	Standard / •	- / •
Accessories & Apps						
Smart Connect	•	•	•	•	•	-
Smart Transmitter	•	•	•	•	•	-
Smart Remote	•	•	•	•	•	-
Smart Key	•	•	•	•	•	•
CROS RIC	•	•	•	-	-	-
Smart Remote App	•	•	•	•	•	•
Smart Connect App	•	•	•	•	•	-
Rexton App	•	•	•	•	•	•





ACCESSORIES & APPS



THE PRACTICAL SMART TRANSMITTER 2,4

THE BLUETOOTH® TRANSMITTER FOR
TELEVISIONS AND STEREO SYSTEMS



The Smart Transmitter 2,4 allows direct transfer of stereo sound from audio sources such as television sets and stereo systems to the hearing systems. Simply place the hearing systems on the transmitter and the pairing process begins. The transmitter can be paired with multiple hearing systems at the same time for shared enjoyment with friends.

SPECIAL FEATURES

- No separate remote control needed
- Connecting/pairing is particularly easy
- Several hearing system users can connect simultaneously
- Analog and digital audio inputs (Cinch and TOSLINK)

THE COMPACT SMART MIC

DIRECT AUDIO STREAMING FROM
PRACTICALLY ANY BLUETOOTH-ENABLED
ANDROID SMARTPHONE



The new Smart Mic enables hands-free telephone calls without your customer needing to pick up his/her smart-phone. In challenging hearing environments it can easily be used as an additional microphone.

SPECIAL FEATURES

- Direct streaming from Android smartphones
- Listen to music in stereo quality
- Can be used as a conference microphone

THE COMFORTABLE CHARGING STATION



The new Charging Station: Not only easy to use for your customers, but also great for the environment. A full charge easily carries through the day, even with streaming music. And - if needed - a quick charge option ensures additional hours of functionality – with only thirty minutes charging time.

SPECIAL FEATURES

- New design with case cover
- Drying function during charging
- Automatic switch-off guarantees overcharge protection
- Inductive charging
- Standard power supply with micro-USB connector

CHARGE ON-THE-GO



Portable charger for StyleLine
and StyleLine BT

M-Core SR hearing aids feature rechargeable lithium- ion batteries for hassle-free power. The included charging case enables on-the-go top-ups, for up to 4 days of use – so your customers can easily enjoy a relaxed weekend trip without needing a power socket.

SPECIAL FEATURES

- Can be charged wirelessly with Qi standard
- Allows additional 3 full charges on the go
- Charging Case recharges fully in just three hours



THE FAST SMART LI-ION POWER

FAST, INDUCTIVE CHARGING



Once the hearing systems are placed in the charging station, the charging process starts automatically. Charging takes place inductively, that is, without direct contact. A complete charge takes approximately 4 hours. Even after just 30 minutes, the hearing systems are ready for another 7 hours of use. The three LEDs display the charging status of the hearing systems.



SPECIAL FEATURES

- Drying function during charging
- Hearing systems automatically switch off when they are placed in the charging station and turn on when they are removed
- Automatic switch-off guarantees overcharge protection
- Standard power supply with micro-USB connector
- Compatible with standard power banks

THE EFFECTIVE SMART POWER

THE ALTERNATIVE TO CONVENTIONAL HEARING SYSTEM BATTERIES



With Smart Power, constant changing hearing system batteries is a thing of the past. Smart Power charges the hearing system batteries reliably and quickly. This avoids the need to fit new batteries constantly.



SPECIAL FEATURES

- Practical charging station
- Environmentally-friendly
- Charges and dries hearing systems overnight

THE COMFORTABLE SMART CONNECT

WIRELESS CONNECTION WITH MANY COMMUNICATION MEDIA



Comfortably listen to music, watch television or make phone calls, also over the Internet – it's all possible. Smart Connect detects the hearing system and pairing takes place automatically. Operation is simple and convenient with the multifunction button or alternatively with the Rexton App.



SPECIAL FEATURES

- Easily activated in Connexx, ready for use without programming
- Simple and quick pairing with the hearing systems
- Simple operation with the multi-function button
- Audio input
- Can also be used hands-free with the integrated microphone

THE VERSATILE SMART TRANSMITTER

CONNECTS NON-BLUETOOTH-ENABLED AUDIO SOURCES WITH SMART CONNECT



As an accessory for Smart Connect, Smart Transmitter connects electronic devices that don't have a Bluetooth® function with the Smart Connect. The push buttons on the top ensure easy and intuitive operation while the LEDs on the front indicate the current status of the transmitter.



SPECIAL FEATURES

- Modern infotainment with Smart Transmitter and Smart Connect
- Directly connected to audio sources
- 3.5 mm jack socket and optical TOSLINK connection

THE EFFECTIVE SMART REMOTE

ATTRACTIVE AND FUNCTIONAL



This remote control can be used with all hearing systems with a wireless function. The timeless design, the large push buttons together with its clearly arranged display are just some of its user-friendly features. Smart Remote offers the following standard functions: changing hearing programs, setting the volume, switching the hearing systems on and off.



SPECIAL FEATURES

- Pre-programmed with standard functions
- Buttons can be customized
- Alarm function
- Service information can be programmed:
Contact information and service appointment reminder
- Adjustable volume of the tinnitus function

THE HANDY SMART KEY

SMALL AND PRACTICAL



With Smart Key, programs and volume can be easily adjusted. A keylock prevents buttons from being inadvertently pressed. A status LED indicates whether the remote control is enabled.



SPECIAL FEATURES

- Small and handy
- Control via high-frequency sounds

THE COMPREHENSIVE REXTON APP

ONE FOR ALL



The Rexton App allows your customers to control their hearing system conveniently. In MOTION CORE and MyCore hearing aids, telephone calls are forwarded directly from the iPhone and take priority over other audio signals; streamed music can be listened to in stereo, directly and without any additional device. With TruCore hearing aids and the Smart Connect you can stream music and telephone calls through this app.

The free Rexton App turns a smartphone into a discreet remote control for the hearing systems. With this app, the hearing systems can be customized and operated for the following settings: hearing program, volume and tone, without the need for any additional accessories.



SPECIAL FEATURES

- Hearing programs and volumes can be adjusted discreetly
- The volume can be set individually for the right and left hearing system
- The Rexton App creates an environment profile for its user. This renders optimization easier for you when it comes to fine-tuning and after sales care

Download the Rexton App now!





RELY ON REXTON

It's a tough world out there, and even tougher when your hearing starts to go. At Rexton, we won't let hearing loss hold you back. Good hearing is critical: From getting the job done, to getting home safely and being there for the people who count on you, we know what's at stake. That's why we work hard to deliver proven hearing technology in the most reliable way. We understand what reliability means in real life, and make practical, easy-to-use products you can count on. Our hearing aids have been getting the job done since 1955, so whatever life throws at you, you can rely on Rexton.

Sivantos GmbH
Henri-Dunant-Straße 100
91058 Erlangen, Germany

www.rexton.com

The Bluetooth® word mark and logos are owned by the Bluetooth SIG, Inc. and any use of such marks by Sivantos GmbH is under license. Other trademarks and trade names are those of their respective owners. Apple and the Apple logo are trademarks of Apple Inc., registered in the U.S. and other countries. App Store is a service mark of Apple Inc. Google Play is a trademark of Google Inc. "Made for iPhone" and "Made for iPad" mean that an electronic accessory has been designed to connect specifically to iPhone or iPad, respectively, and has been certified by the developer to meet Apple performance standards. Apple is not responsible for the operation of this device or its compliance with safety and regulatory standards. Please note that the use of this accessory with iPhone or iPad may affect wireless performance. iPad, iPhone are trademarks of Apple Inc., registered in the U.S. and other countries.

Order No. 04140-99T1-7600 © 2019 Sivantos GmbH