ordering information

Low profile Leo+ baby offers the possibility to treat small arteries down to \varnothing 1,5 mm

	Reference	Vessel Ø	Stent's length	Unconstro	nined stent	Compatible
		(mm)	at nominal Ø (mm)	Ø (mm)	Length (mm)	delivery catheters
Leo+ baby 2,0	LEO. 2,0x12		12	2,8	8	
	LEO. 2,0x18	1,5 to 2,5	18		12	VASCO+10
	LEO. 2,0x25*		25		16	
Leo+ baby 2,5	LEO. 2,5x12		12	3,2	8	
	LEO. 2,5x18		18		12	
	LEO. 2,5x25	2,0 to 3,1	25		16	VASCO+10
	LEO. 2,5x30*		29		19	
	LEO. 2,5x35*		35		23	
Leo+ baby 3,0	LEO.3.0x12	2,5 to 3,6	12	2.0	9	
	LEO. 3,0x18		18		13	VASCO+10
	LEO. 3,0x25	2,3 10 3,0	25	3,8	16	VA3CO+10
	LEO. 3,0x35		35		24	
Leo+ 3,5	LEO. 3,5x12	_	12	4,4	9	
	LEO. 3,5x18		18		12	
	LEO. 3,5x25	3,10 to 4,25	25		17	VASCO+21
	LEO. 3,5x30	3,10104,23	30		19	VAJCOTZI
	LEO. 3,5x35		35		24	
	LEO. 3,5x50		50		36	
Leo+ 4,5	LEO. 4,5x15		15		12	
	LEO. 4,5x20		21		15	
	LEO. 4,5x25		26		18	
	LEO. 4,5x30	4,25 to 5,30	30	5,4	21	VASCO+25
	LEO. 4,5x40		40		28	
	LEO. 4,5x50		51		37	
	LEO. 4,5x75		75		55	
Leo+ 5,5	LEO. 5,5x25		26	6,7	18	
	LEO. 5,5x30		31		21	
	LEO. 5,5x35	5,30 to 6,50	36		26	VASCO+28
	LEO. 5,5x50	3,00 10 0,00	52		35	W 10COTZ0
	LEO. 5,5x60		60		42	
	LEO. 5,5x75		75		51	

*Manufactured on demand

Reference: 1. Internal data

The self-expandable LEO+/LEO+ Baby stent is designed for the treatment of intracranial aneutysms in association with embolization coils. Class III CE0297 in compliance with Medical Device Directive (MDD 93/42/EEC amended by 2007/47/EC). Manufactured by BALT Extrusion SAS. Carefully read the instructions for use before use. First CE marking:2007 (LEO+),2012 (LEO+Baby). VASCO+ is a reinforced micro-catheter intended for injection of diagnostic or therapeutic product, to position pushable coils "SPIRALES" or detachable coils especially the ones of MDS * mechanical detachment system *, for the use of the self-expanding stent LEO+ or SILK+. Class III CE0297 in compliance with Medical Device Directive (MDD 93/42/EEC amended by 2007/47/EC). Manufactured by BALT Extrusion SAS. Carefully read the instructions for use before use. First CE marking:2004.

The content of this document, in particular data, information, trademarks and logos is BALT SAS and affiliates' sole property. © 2019 BALT SAS and affiliates, all rights reserved. All representation and/or reproduction, whether in part or in full, is forbidden and would be considered a violation of BALT SAS and affiliates' copyrights and other intellectual proprietary rights. This document with associated pictures is non-contractual and is solely dedicated to healthcare professionals and BALT's distributors (BALT's supplier's distributors). It cannot be distributed or given to patients. The products commercialized by BALT shall exclusively be used in accordance with the instructions for use included in the boxes.

DC015GB (01/2019)

Bal

10, rue de la Croix Vigneron, 95160 Montmorency France

Tel.: +33 (0)1 39 89 46 41 Fax: +33 (0)1 34 17 03 46

www.balt.fr





eo+& leo+ baby

Designed for the treatment of intracranial aneurysms in association with embolization coils

More than

40 000

Leos deployed around the world since 2003





Flexibility thanks to the nickel titanium materials

Conformability ven by the sliding-cell technology

Smooth wall apposition as a result of the rounded short flared ends

Radial force for optimal coil mass support

Intraluminal support & neck coverage thanks to the 16 braided wires

Accurate positioning

Visibility

two helical markers on the entire body of the stent

Resheathability

up to 90% of deployed length


Flow diversion effect

Bouillot et al., Computational fluid dynamics with stents: quantitative comparison with particle image velocimetry for three commercial off the shelf intracranial stents J NeuroIntervent Surg 2016; 8:309-315"





braided wires

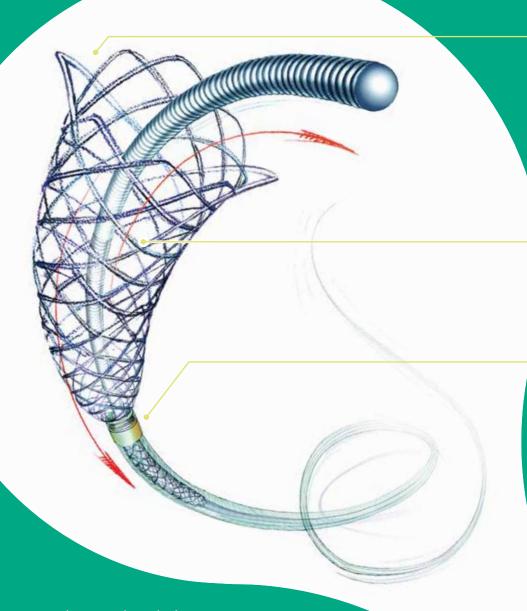

resheathable
up to %
of deployed length¹

Range details – leo+ baby Technical features

		Deference	Vessel Ø		Unconstr	D. P		
Ke		References		Stent's length at nominal Ø (mm)	Ø (mm)	Length (mm)	Delivery catheter	
LEO+ BABY 2,0		LEO2,0X12		12		8		
	LEO2,0X18	1,5 to 2,5	18	2,8	12	VASCO+10		
	LEO2,0X25*		25		16			
LEO+ BABY 2,5	LEO2,5X12	2,0 to 3,1	12	3,2	8			
	LEO2,5X18		18		12	VASCO+10		
	LEO2,5X25		25		16			
		LEO2,5X30*		30		19		
	LEO2,5X35*		35		23			
LEO+ BABY 3,0		LEO3,0X12	2,5 to 3,6	12	2.0	9	VASCO : 10	
		LEO3,0X18		18		13		
	7	LEO3,0X25		25	3,8	16	VASCO+10	
		LEO3,0X35		35		24		

LEO+/LEO+ Baby, self-expandable intracranial stent intended for the treatment of intracranial aneurysms in association with embolization coils. Class III CE0297 in compliance with Medical Device Directive (MDD 93/42/EEC amended by 2007/47/EC). Manufactured by BALT EXTRUSION SAS. Carefully read the instructions for use before use. First CE marking 2007 (LEO+) 2012 (LEO+Baby). The content of this document, in particular data, information, trademarks and logos is BALT SAS and affiliates' sole property. ©2022 BALT and affiliates all rights reserved. Consequently, all representation and/or reproduction, whether in part or in full, is forbidden and would be considered a violation of BALT SAS and affiliates' copyrights and other intellectual proprietary rights. This document with associated pictures is non-contractual and is solely dedicated to BALT's distributors (BALT's supplier's distributors). The products commercialized by BALT shall exclusively be used in accordance with the package inserts which have been updated and included in the boxes. The distribution, the sale and the use of BALT's products is under distributor's liability.



INTERNAL USE -CONFIDENTIAL

Range details – leo+

nical features	REFERENCE	VESSELØ	STENT's LENGTH AT NOMINAL Ø (mm)	UNCONSTRAINED STENT		
	REFERENCE	VESSELØ		Ø (mm)	LENGTH (mm)	DELIVERY CATHETER
	LEO3,5X12	3,10 to 4,25	12	4,4	9	VASCO+21
	LEO3,5X18		18		12	
	LEO3,5X25		25		17	
LEO+ 3,5 -	LEO3,5X30		30		19	
	LEO3,5X35		35		24	
	LEO3,5X50		50		36	
	LEO4,0X15	3,75 to 4,75	17	4,0	12	VASCO+21
	LEO4,0X20		21		15	
LEO+ 4,0	LEO4,0X25		26		18	
	LEO4,0X30		30		21	
New	LEO4,0X35		36		25	
ffering a bigger diameter	LEO4,0X40		41		29	
till deliverable through a	LEO4,0X50		52		37	
Vasco+21	LEO4,5X15		12		12	
	LEO4,5X20	4,25 to 5,30	21	5,4	15	VASCO+25
	LEO4,5X25		26		18	
LEO+4,5 -	LEO4,5X30		30		21	
	LEO4,5X40		40		28	
	LEO4,5X50		51		37	
	LEO4,5X75		75		55	
	LEO5,5X25	5,30 to 6,50	26	6,7	18	VASCO+28
LEO+ 5,5 -	LEO5,5X30		31		21	
	LEO5,5X35		36		26	
	LEO5,5X50		52		35	
	LEO5,5X60		60		42	
	LEO5,5X75		75		51	



INTERNAL USE -CONFIDENTIAL