

Ordering Information

Tests for choices	Specimen	Description	Catalog No.	Format
MOP、MET、THC、MDMA、KET、K2、AMP、BAR、BUP、BZO、COC、COT、EDDP、FTY、MQL、MTD、OPI、OXY、PCP、PPX、TCA、TRA、CAF、ETG、MCAT	Urine	Multi-Drug Panel	DOA-A*4a	2 -12 drugs

⚠ Note: For drug combination, please contact us for more details.

Drug of Abuse | Urine-CUP

The Pocguide™ Multi-Drug Test Cup was designed to facilitate urine collection and drug testing in one cup. This fully functioning cup design eliminates the handling of urine. The cup is a popular choice for convenient drug testing in workplace and home settings.

- **Flexibility**
 - 2 to 18 drugs per cup;
 - Custom configurations available;
 - Wide choice of test combination: 25 drugs assays to choose from;

▪ Convenient Test Procedure


- No specimen handing required
- Eliminates test reagent contamination
- Convenient visual interpretation
- On-site preliminary screening

▪ Performance


- Accurate results:high correlation with GC/MS
- Fast result:5 minutes test time

▪ How to use CUP

- Test must be in room temperature (15°C to 30°C).
- After the urine has been collected, re-cap the cup and place the test cup on a flat surface.
- Peel the label from right to left and read the result within 5 minutes. Do not read results after 10 minutes.

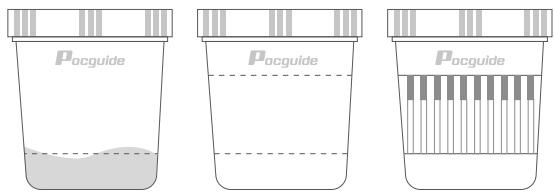

minimum urine level remove this label to view results

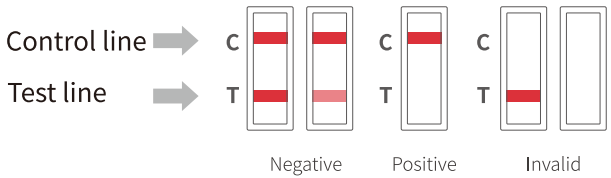

Ordering Information

Tests for choices	Specimen	Description	Catalog No.	Format
MOP、MET、THC、MDMA、KET、K2、AMP、BAR、BUP、BZO、COC、COT、EDDP、FTY、MQL、MTD、OPI、OXY、PCP、PPX、TCA、TRA、CAF、ETG、MCAT	Urine	Multi-Drug CUP	DOA-A*7a	2 -18 drugs

⚠ Note: For drug combination, please contact us for more details.

Drug of Abuse | Saliva Test

DOA Oral Fluid Rapid Tests

▪ Performance

- Accurate results
- Results in 10 minutes
- High sensitivity and specificity
- Tested for cross-reactivity against most common drugs and other substances

▪ Advantage

- Non-invasive
- No bathroom needed
- Easy and rapid sample collection
- No specimen dilution or adulteration

Ordering Information

Tests for choices	Sensitivity	Specimen	Description	Catalog No.	Format
MOP/OPI	10ng/ml	Saliva	Multi-Drug Saliva	DOA-J*5a	2 -12 drugs
MET	25ng/ml				
THC	20ng/ml				
AMP	25ng/ml				
BAR	25ng/ml				
BUP	10ng/ml				
BZO	5-10ng/ml				
COC	20ng/ml				
FTY	10ng/ml				
MTD	30ng/ml				
OXY	40ng/ml				
TRA	25ng/ml				

⚠ Note: For drug combination, please contact us for more details.



Email:international@aichek.com Tel:+8657186963020
Add:3rd Floor,Building 9,Hexiang Technology Center, Xiasha Street, Qiantang District,
310018 Hangzhou,Zhejiang,P.R.China



Drug of Abuse
Rapid Diagnostic Solutions



The Pocguide™ Drugs of Abuse Rapid Tests portfolio provides fast and accurate drug screening testing options on urine and oral fluids in a variety of convenient test device formats. Our rapid tests are used in a multitude of settings such as hospitals, clinics, physician offices, criminal justice systems, home testing, education, immigration systems and workplaces such as mining, construction, transport and aviation.


Product Overview

The Pocguide™ DOA Rapid Tests portfolio offer Single and Multi-Drug Test Devices, Drug Tests with Integrated Cup Systems,Panel Tests and Saliva Tests. Our rapid test technology is based on lateral flow chromatographic immunoassay for the simultaneous, qualitative detection of multiple drugs and metabolites in human urine or saliva at selected cut-off values.

Drug of Abuse | Urine-Device

DOA Single-drug device Tests:A range of 24 common drugs are now available with a single device to facilitate single drug testing.

Test Procedure:

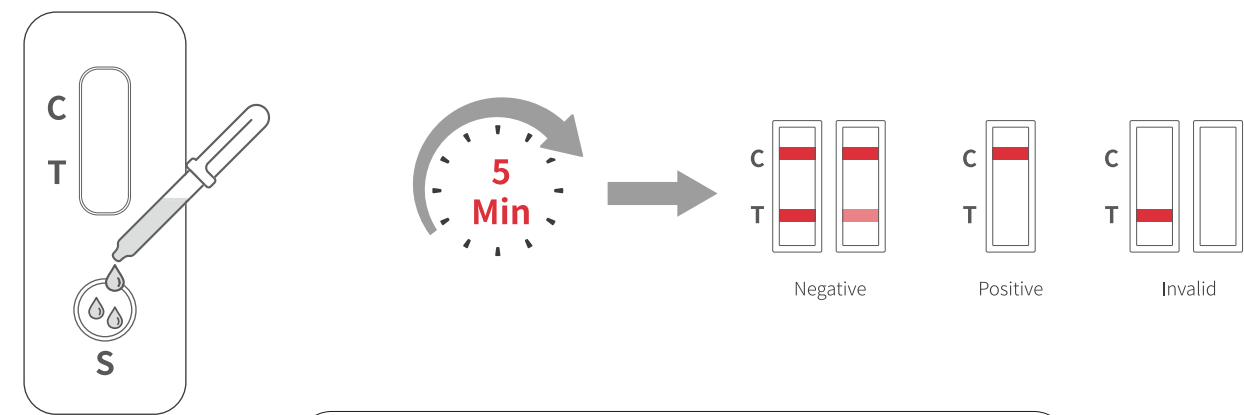

- DEVICE TEST
- Step1.

Add 3 full drops of urine
- Step2.

Read Result at 5 minutes

Ordering Information

Product Description	Specimen	Format	Catalog No.	Sensitivity	Kit Size
MOP One Step Morphine Test Device	Urine	Device	MOP-A02a	300 ng/mL	25
MET One Step Methamphetamine Test Device	Urine	Device	MET-A02a	1000 ng/mL	25
THC One Step Marijuana Test Device	Urine	Device	THC-A02a	50 ng/mL	25
MDMA One Step Ecstasy Test Device	Urine	Device	MDMA-A02a	500 ng/mL	25
KET One Step Ketamine Test Device	Urine	Device	KET-A02a	1000 ng/mL	25
K2 One Step Synthetic Cannabis Test Device	Urine	Device	K2-A02a	25 ng/mL	25
AMP One Step Amphetamine Test Device	Urine	Device	AMP-A02a	1000 ng/mL	25
BAR One Step Barbiturates Test Device	Urine	Device	BAR-A02a	100 ng/mL	25
BUP One Step Buprenorphine Test Device	Urine	Device	BUP-A02a	10 ng/mL	25
BZO One Step Benzodiazepine Test Device	Urine	Device	BZO-A02a	300 ng/mL	25
COC One Step Cocaine Test Device	Urine	Device	COC-A02a	300 ng/mL	25
COT One Step Cotinine Test Device	Urine	Device	COT-A02a	100 ng/mL	25
One Step EDDP Test Device	Urine	Device	EDDP-A02a	100 ng/mL	25
FTY One Step Fentanyl Test Device	Urine	Device	FTY-A02a	20 ng/mL	25
MQL One Step Methaqualone Test Device	Urine	Device	MQL-A02a	300 ng/mL	25
MTD One Step Methadone Test Device	Urine	Device	MTD-A02a	300 ng/mL	25
OPI One Step Opiates Test Device	Urine	Device	OPI-A02a	2000 ng/mL	25
OXY One Step Oxycodone Test Device	Urine	Device	OXY-A02a	100 ng/mL	25
PCP One Step Phencyclidine Test Device	Urine	Device	PCP-A02a	25 ng/mL	25
PPX One Step Propoxyphene Test Device	Urine	Device	PPX-A02a	300 ng/mL	25
TCA One Step Tricyclic Antidepressants Test Device	Urine	Device	TCA-A02a	1000 ng/mL	25
TRA One Step Tramadol Test Device	Urine	Device	TRA-A02a	100 ng/mL	25
CAF One Step Caffeine Test Device	Urine	Device	CAF-A02a	1000 ng/mL	25
ETG One Step Ethyl Glucuronide Test Device	Urine	Device	ETG-A02a	500 ng/mL	25
MOP/MET/KET/MDMA/THC One Step Test	Urine	Device	DOA-A155a	See Insert	25

Drug of Abuse | Urine-Panel

Our Multi-Drug Test Panel is a test for the simultaneous detection of 2-12 drugs of abuse in a single urine specimen.

Features and Benefits

- Easy to use: just dip then read;
- Results ready in 5 minutes;
- Accurate results: high correlation with GC/MS


Flexibility

- 2 to 12 drugs per panel;
- Custom configurations available;
- Wide choice of test combination: 25 drugs assays to choose from;



Test Procedure

- Allow the test panel and urine specimen to room temperature (15-30°C) prior to testing.
- Open the sealed pouch by tearing along the notch. Remove the test device from the pouch.
- Hold the one side of the device with one hand. Use the other hand to pull out the cap and expose the absorbent end.
- Immerse the absorbent end into the urine sample about 10 seconds. Make sure that the urine level is not above the “MAX” line printed on the front of the device.
- Lay the device flat on a clean, dry, non-absorbent surface.
- Read the result at 5 minutes. Do not read after 10 minutes.

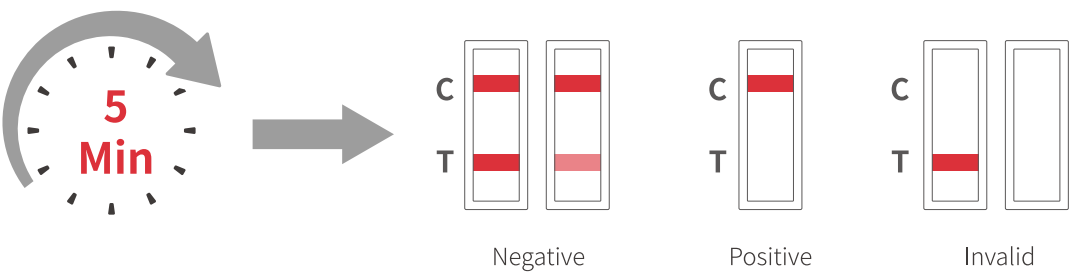

Note:For detailed instruction, limitation and warning, please refer to product insert.