

Microcure (Suzhou) Medical Technology Co.,Ltd.

Web: https://www.microcuremed.com/

Add: 1-101,188 Jialingjiang Road, Suzhou New District, Jiangsu, P.R.China.

Contact: Jerrish Young (Business Development Manager)

TEL: +86-0512-62916116

Mobile: +86-13681881006

Email: export@microcure.com.cn



LinkedIn



Scan code to visit our website



INS

microcure

IN PERFECTION
WE BELIEVE





Company profile	/01	Disposable linear stapler	/19
Corporate culture	/03	Disposable trocar	/21
Disposable endoscopic linear cutter		Disposable Trocars For Thyroid Surgery	/25
stapler and reloads	/05	Veress Needle	/26
Disposable skin stapler	/09	Disposable incision protection sleeve	/27
Disposable circular stapler	/13	Surgery retractor system	/31
Disposable hemorrhoid stapler	/15	Disposable sample retrieval bag	/35
Disposable reloadable linear cutter stapler and reloads	/17	Disposable ligation loop	/37
		Certificate	/39



MICROCURE

We are committed to pursuing excellence in both design and provision of minimally invasive surgical instruments and patient-centered services.







sion

An enterprise to be accomplished by people of empathy, who are striving to make the world healthier

Credo In perfection we believe

Strategy
Catch up via differentiation


DISPOSABLE ENDOSCOPIC LINEAR CUTTER STAPLER AND RELOADS

Applications

Applicable for resection, transection and anastomosis of lung, bronchial tissue, stomach and intestine in open or endoscopic surgery.





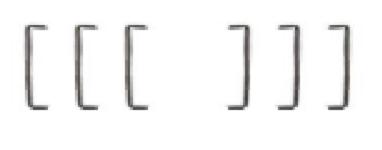


B-shape Staple

Product Advantages

01 Curved tip

Not only provides better visualization to ensure the safety of resection and anastomosis of blood vessels, but also allows blunt dissection of tissues or organs.

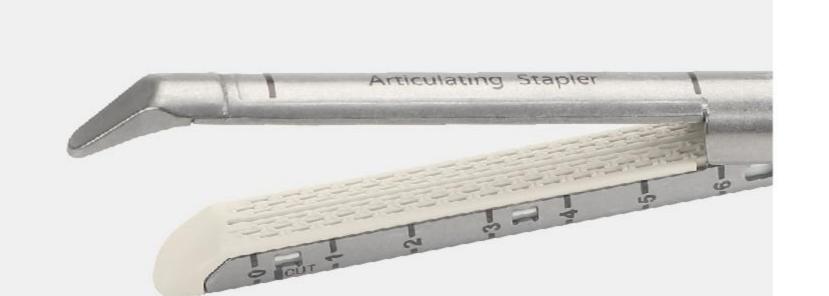



02 "Droplets" nail groove

Optimize B-shape staple formation and delivers superior hemostasis and pressure resistance.

03 Unibody anvil

Render bite force adequate and complete compression.





04 Three-row unequal height stitching design

Progressive clamping of the tissue to ensure staple formation, facilitate lateral drainage of tissue fluid, easier compression to the right height, ensure staple formation, provide wider adaptability of tissue thickness, tissue healing, better closure and haemostasis.



05 Single-hand operation

Make reloads to easily clamp and release, free up surgeons' hands.



06 Rotating system

The marble-type rotary nob delivers more smooth, stable and easy rotation, and the large wave-shaped knob for easy one-handed operation.

07 Different specifications

The gun has 3 models total, and is compatible with reloads of all specifications, which also include different models according to their colors and heights, helping to deliver superior performance according to a range of thickness.




Order Information

STAPLER	ORDER CODE	LENGTH(mm)	SHAFT DIAMETER(mm)
	WZULS-A-S	60	12
	WZULS-A-L	160	12
STAPLER	WZULS-A-XL	260	12
BODY	WZULS-B-S	60	12
	WZULS-B-L	160	12
	WZULS-B-XL	260	12
RELOADS	MODELS	SPECIFICATIONS	COLOR OF RELOADS
	WZULS-P1	30-2.5	white
	WZULS-P1	45-2.5	white
	WZULS-P1	45-3.5	blue
STRAIGHT RELOADS	WZULS-P1	45-4.8	green
KELOADS	WZULS-P1	60-2.5	white
	WZULS-P1	60-3.5	blue
	WZULS-P1	60-4.8	green
	WZULS-P2	30-2.5	white
	WZULS-P2	45-2.5	white
	WZULS-P2	45-3.5	blue
ARTICULATING	WZULS-P2	45-4.8	green
RELOADS	WZULS-P2	60-2.5	white
	WZULS-P2	60-3.5	blue
	WZULS-P2	60-4.8	green
	WZULS-P3-30B	2.0-2.5-3.0	brown
	WZULS-P3-30BE	2.0-2.5-3.0	brown
	WZULS-P3-30P	3.0-3.5-4.0	purple
	WZULS-P3-45B	2.0-2.5-3.0	brown
TRI STAPLE	WZULS-P3-45BE	2.0-2.5-3.0	brown
	WZULS-P3-45P	3.0-3.5-4.0	purple 🔲
	WZULSP3-45H	4.0-4.5-5.0	black
	WZULS-P3-60B	2.0-2.5-3.0	brown
	WZULS-P3-60P	3.0-3.5-4.0	purple 🔲
	WZULS-P3-60H	4.0-4.5-5.0	black

*WZULS-A is with rubber inlay and WZULS-B is without rubber inlay.
*OEM services are available.

microcure 108 07 microcure

DISPOSABLE SKIN STAPLER

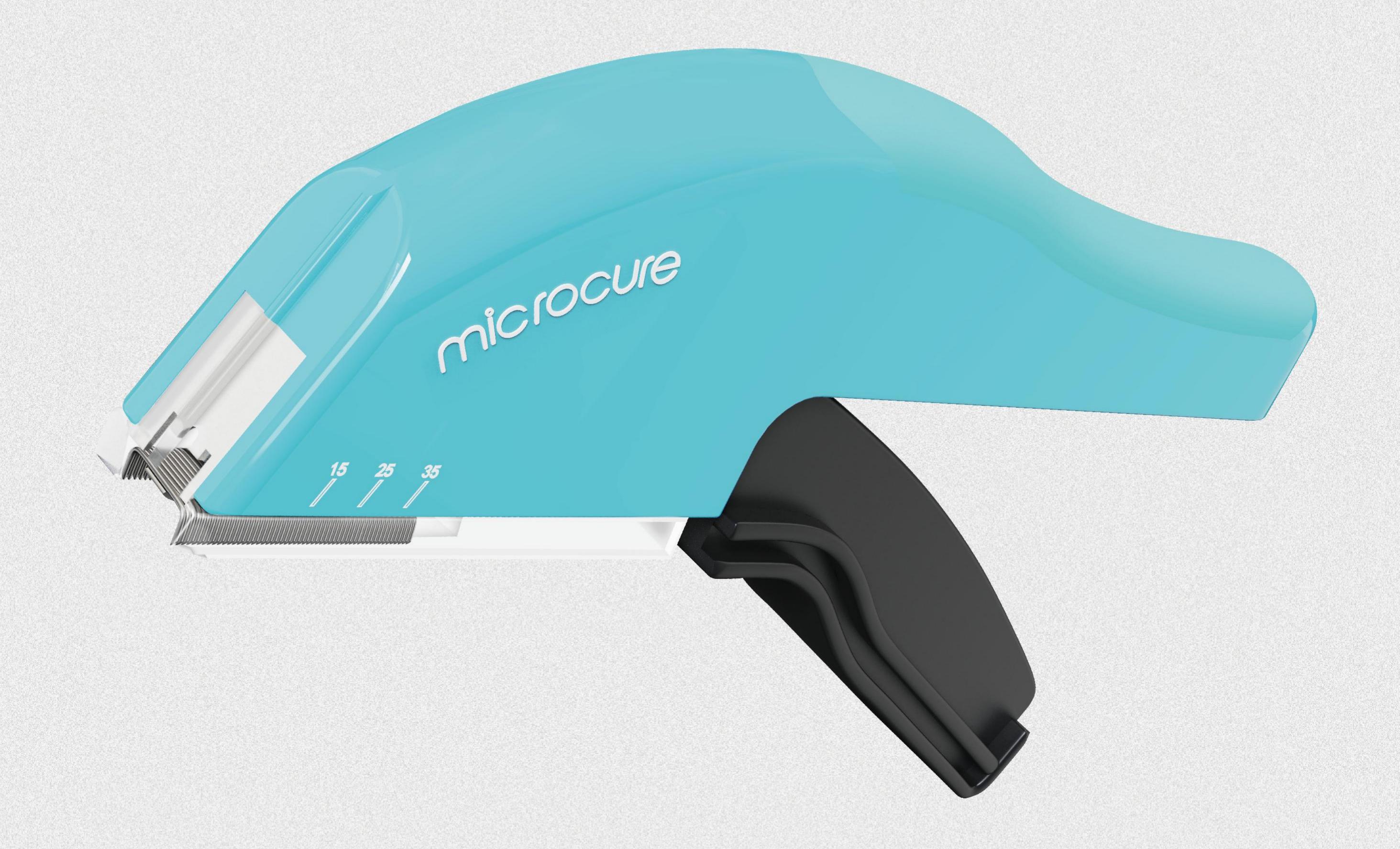

Applications

Applicable to the closure of surgical and traumatic wounds (straight and roughly straight shapes)

APPLICABLE DEPARTMENT

01 EMERGENCY DEPARTMENT

02 GENERAL SURGERY

03 ANIMAL SURGERY

04 CARDIOTHORACIC SURGERY

Product Advantages

01 Simple to use

Provide simple operation, lightening the burden on surgeons.

35W

02 Fast suture

Eases patients' pains and improves operation efficiency.

nicrocure i i i i

03 Quick recovery

Minimizes the incision scar and postoperative recovery time.




Stapler

01 Evert and approach skin incisions.


02 Release trigger to achieve the suturing.

03 Place stapler on the center of everted skin incisions and squeeze the trigger completely.

Remover

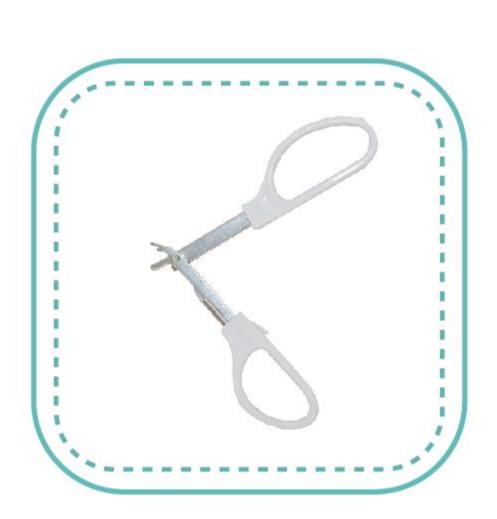


01 Insert the front part of the remover under suture nails, and slide the nails into the notch in the lower jaw.

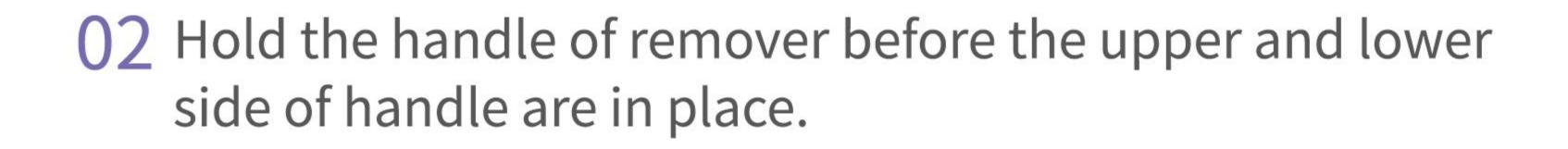




03 Confirm the suture nails deformed and separeted from the skin, and remove the nails.


Product Demo

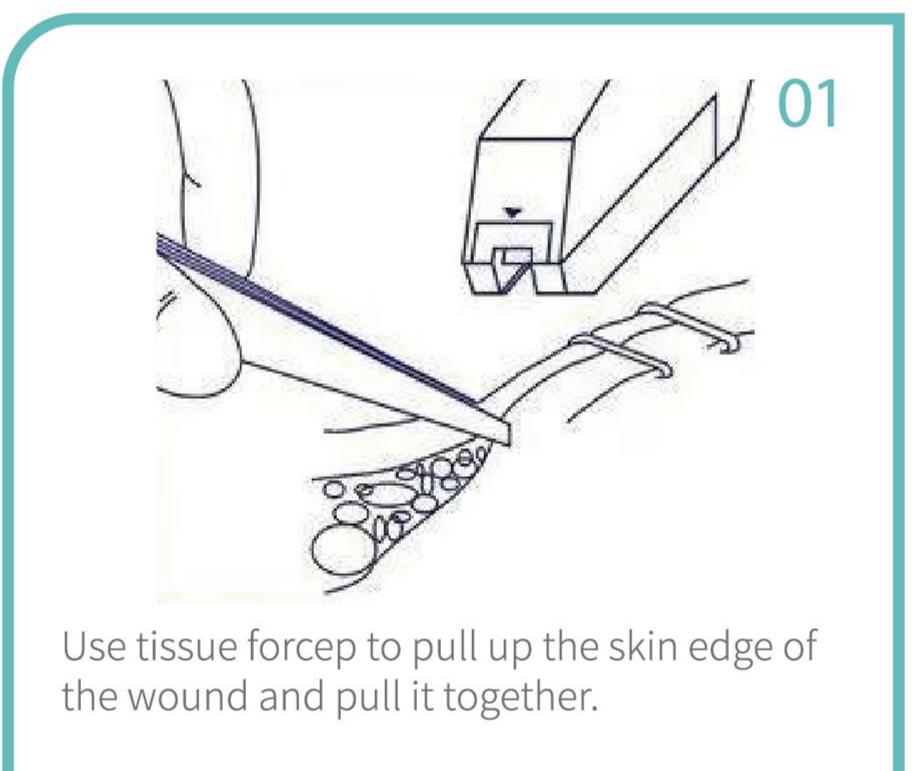

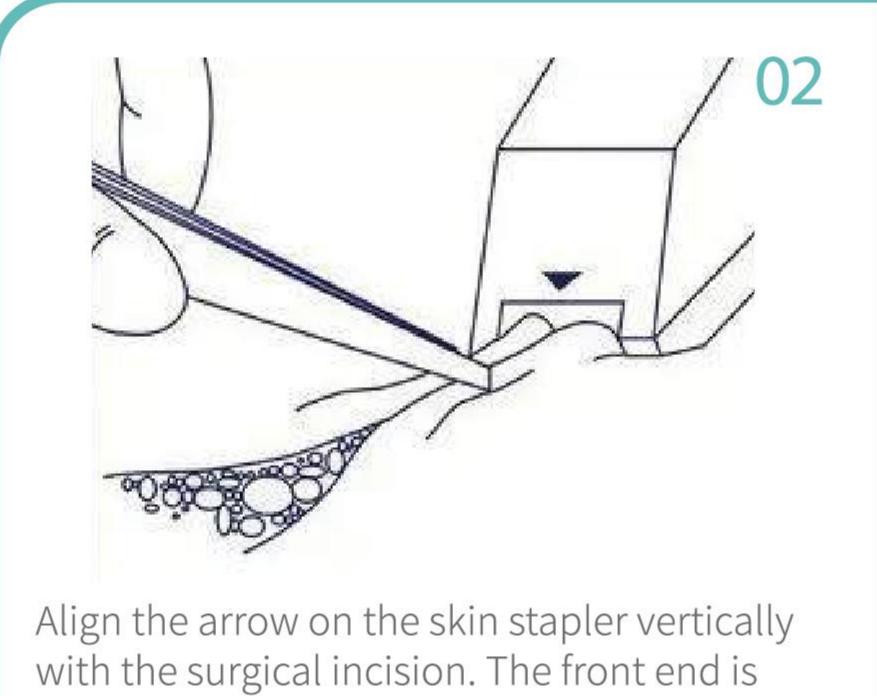

handle is pressed in place.

close to the skin, hold the upper and lower handle tightly, and apply force evenly until the

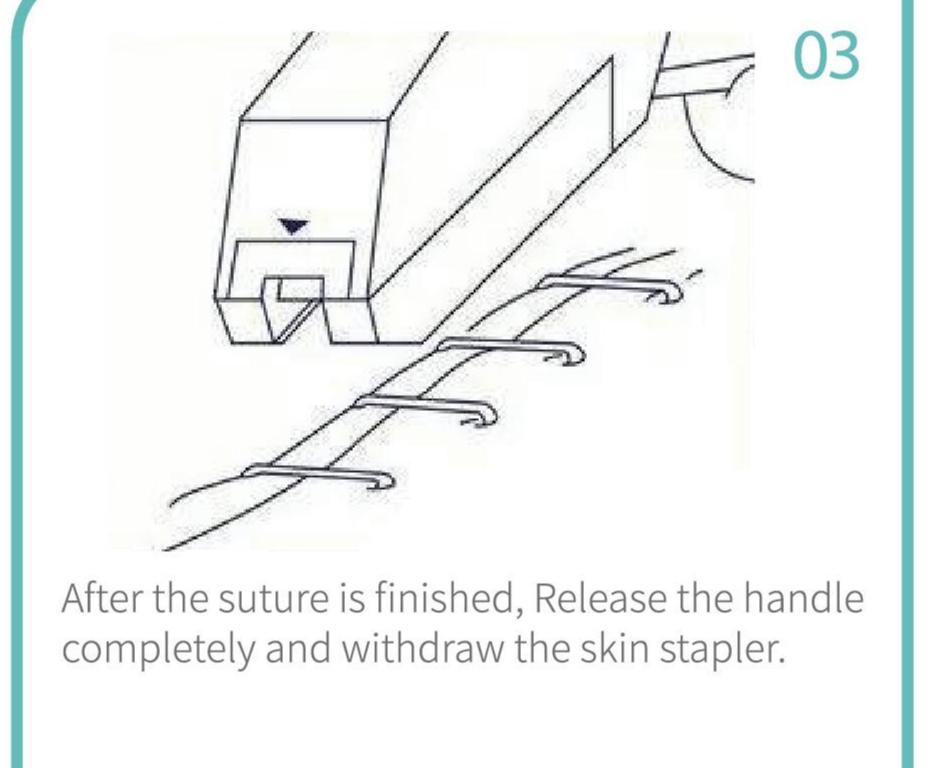

Order Information

MODELS	STAPLER WIDTH (mm)	STAPLER DIAMETER (mm)
WZDSS-A-15N	12	0.6
WZDSS-A-15W	13.8	0.6
WZDSS-A-25N	12	0.6
WZDSS-A-25W	13.8	0.6
WZDSS-A-35N	12	0.6
WZDSS-A-35W	13.8	0.6
WZDSS-A-45N	12	0.6
WZDSS-A-45W	13.8	0.6


NAME	MODELS	
Skin stapler remover	WZDSS-P1	

*OEM services are available.

microcure 12 11 microcure

DISPOSABLE CIRCULAR STAPLER



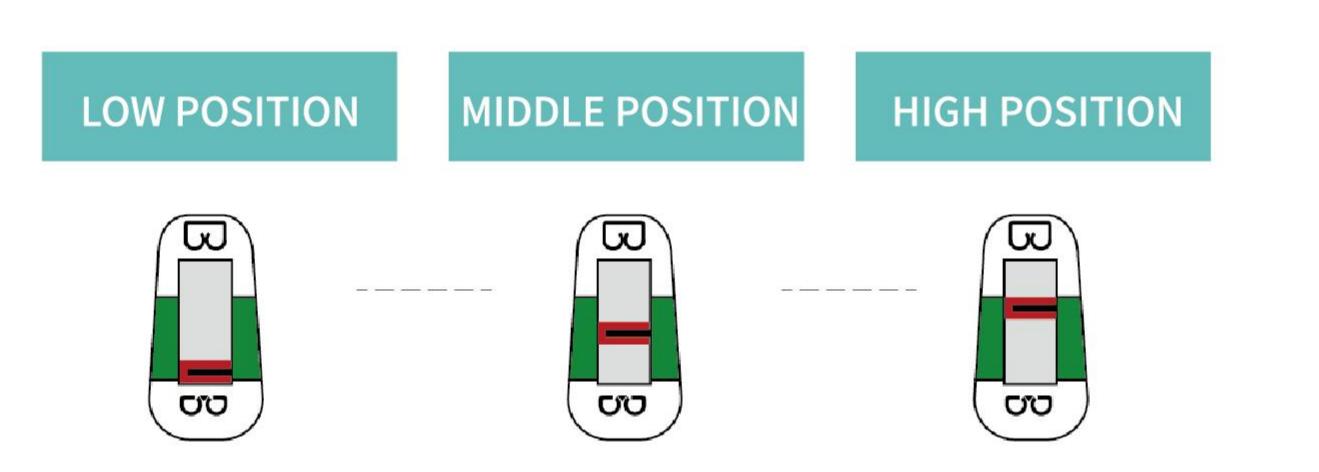

Applications

Applicable to end to end, end to side, side to side anastomosis in digestive tract reconstructive surgery of esophagus, stomach and intestines.

Product Advantages

O1 Minimum suture gap Prevent tissue been crushed.	O2 Four-axis positioning nail Ensure precise staple formation.
O3 Material Titanium staple provides tensile strength to ensure steady anastomosis after surgery.	04 Ring-shape blade Made through special process and cut sharply.
05 Indicator window Provide real-time position of nail.	06 B-shape staple Reduce the risk of postoperative bleeding.

Order Information

MODELS	TOTAL LENGTH (mm)	OUTER DIAMETER (mm)	STAPLE HEIGHT (mm)	STAPLE DIAMETER (mm)
WZCS-A-24	450	24	5	0.28
WZCS-A-26	450	26	5	0.28
WZCS-A-29	450	29	5.2	0.28
WZCS-A-32	450	32.5	5	0.28

^{*}OEM services are available.

^{*}WZCS-A-15 and WZCS-A-33 are available.

DISPOSABLE HEMORRHOID STAPLER

Applications

Applicable for selective resection and suturing of dentate line mucos.



LOW POSITION MIDDLE POSITION HIGH POSITION

VARIOUS ACCESSORIES

VARIOUS ACCESSORIES







Product Advantages

01	Minimum suture gap	
Prev	ent tissue been crushed.	

O2 Ring-shape blade

Made through special process and cut sharply.

O3 Indicator window
Provide real-time position of nail.

O4 Larger cartridge
Provide large scope for more tissue.

05 B-shape staple Reduce the risk of postoperative bleeding.

O6 Reloads
A broad range of specifications meet various surgical requirements.

Order Information

MODELS	TOTAL LENGTH (mm)	OUTER DIAMETER (mm)	STAPLE HEIGHT (mm)	STAPLE DIAMETER (mm)
WZHS-A-32	395	32	3.8	0.28
WZHS-A-33	395	33.5	3.8	0.28
WZHS-A-34	395	34	3.8	0.28


^{*}OEM services are available.

^{*}Stitch piece for anorectal examination can be sold separately.

DISPOSABLE RELOADABLE LINEAR CUTTER STAPLER AND RELOADS

Applications

To use in creation of anastomotic site or closure of incision in recreation surgeries of digestive tract or other resection procedures.

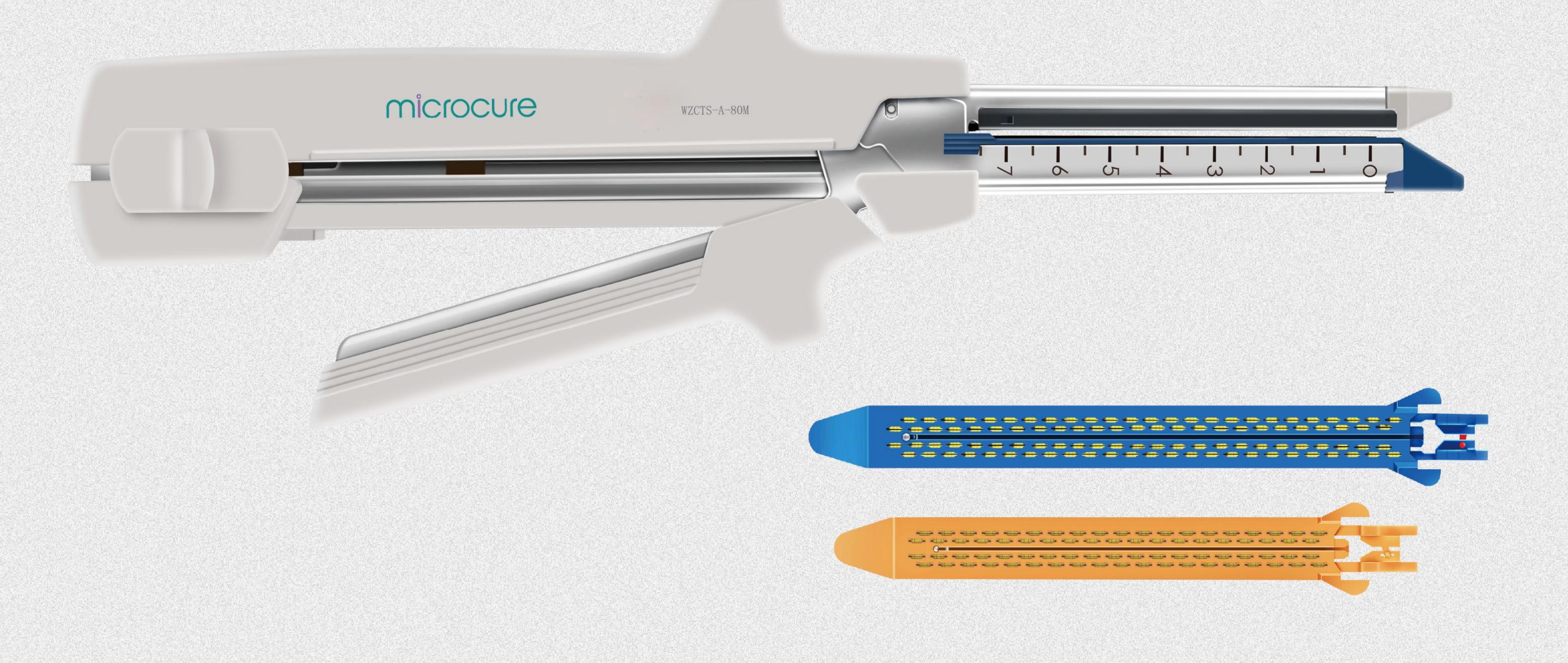

Applicable departments

General surgery, thoracic surgery, hepatobiliary surgery.

Product Advantages

01 Anvil

Anvil with high strength meets more challenging surgical requirements.

02 Length

The length of stapled line is 1.5 times of that of staple, preventing the risk of bleeding.

03 Convenience

Combination between manual and automatic operation along with corresponding sounds.

04 Reloadable

Reloadable design allows 5 replacements of reloads during one surgery, thus reducing costs.

Order Information

Models	Staple diameter (mm)	Staple height (mm)	Staple amount
WZCTS-A-60M	0.2	3.8	60
WZCTS-A-60H	0.2	4.2	60
WZCTS-A-80M	0.2	3.8	80
WZCTS-A-80H	0.2	4.2	80
WZCTS-A-100M	0.2	3.8	100
WZCTS-A-100H	0.2	4.2	100
WZCTS-P1-60M	0.2	3.8	60
WZCTS-P1-60H	0.2	4.2	60
WZCTS-P1-80M	0.2	3.8	80
WZCTS-P1-80H	0.2	4.2	80
WZCTS-P1-100M	0.2	3.8	100
WZCTS-P1-100H	0.2	4.2	100

*OEM services are available. 18

DISPOSABLE LINEAR STAPLER

Applications

To suture incision or open end in digestive tract recreation surgery or viscera resection.





Applicable departments

Applicable for gastrointestinal surgery.

Product Advantages

- 01 B-shape StapleReduce the risk of postoperative bleeding.
- O2 Narrower design
 Reduce the risk of causing damage to peripheral tissue in surgery.
- O4 Reloads with a broad range of models meet various surgical requirements.
- 03 Reset buttonMake operation easy and efficient.

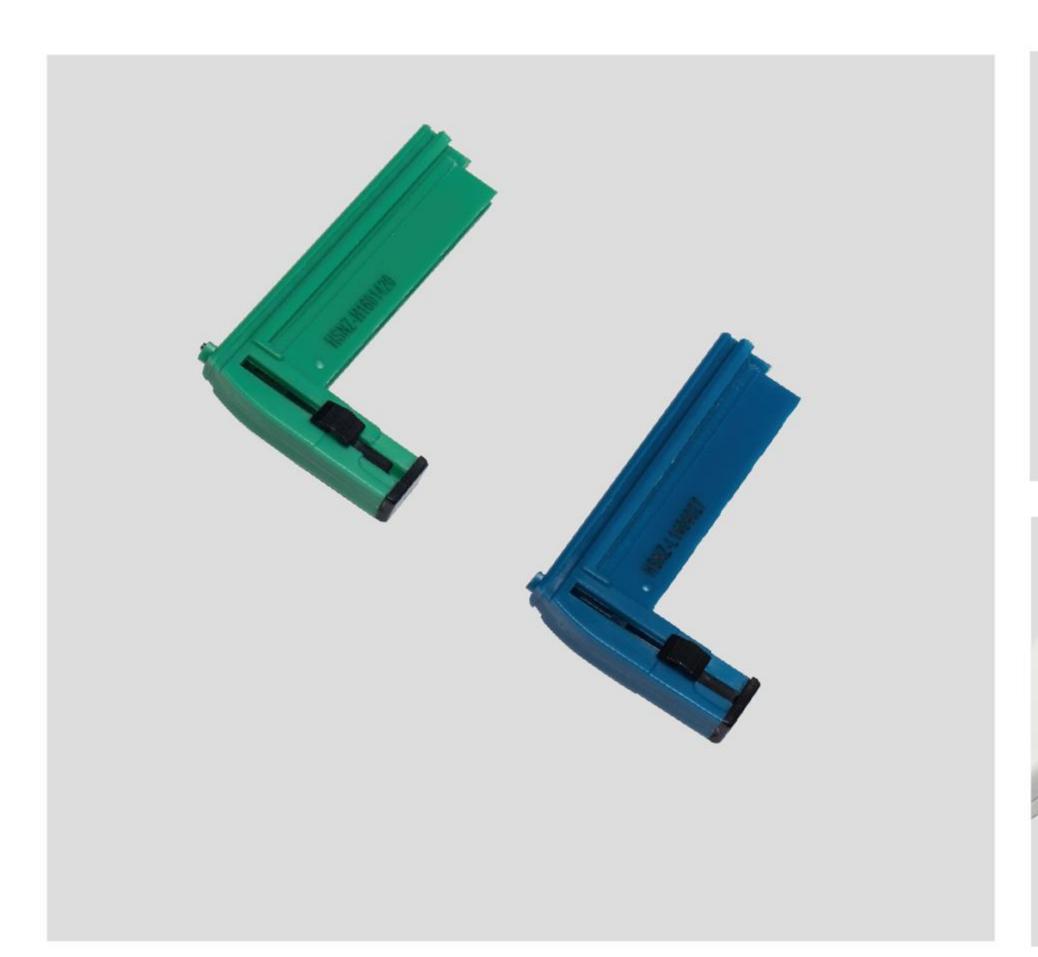

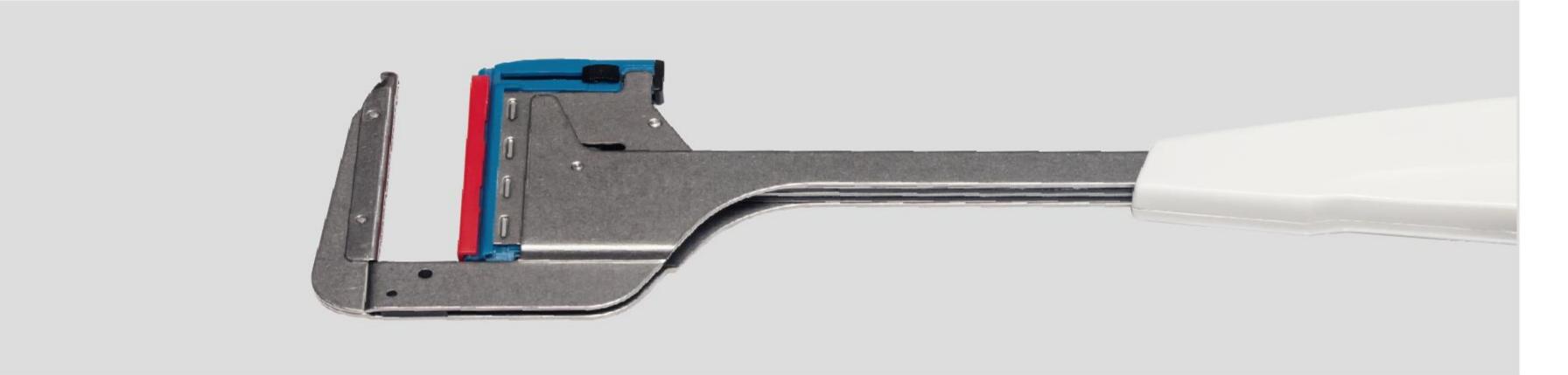






Order Information

MODELS	STAPLE DIAMETER(mm)	STAPLE HEIGHT(mm)	STAPLE AMOUNT	COLOR OF RELOADS
WZLS-A-30L	0.28	3.8	13	blue
WZLS-A-30H	0.28	4.5	13	green
WZLS-A-45L	0.28	3.8	19	blue =
WZLS-A-45H	0.28	4.5	19	green
WZLS-A-60L	0.28	3.8	25	blue
WZLS-A-60H	0.28	4.5	25	green
WZLS-A-90L	0.28	3.8	39	blue =
WZLS-A-90H	0.28	4.5	39	green


*OEM services are available.

DISPOSABLE TROCAR

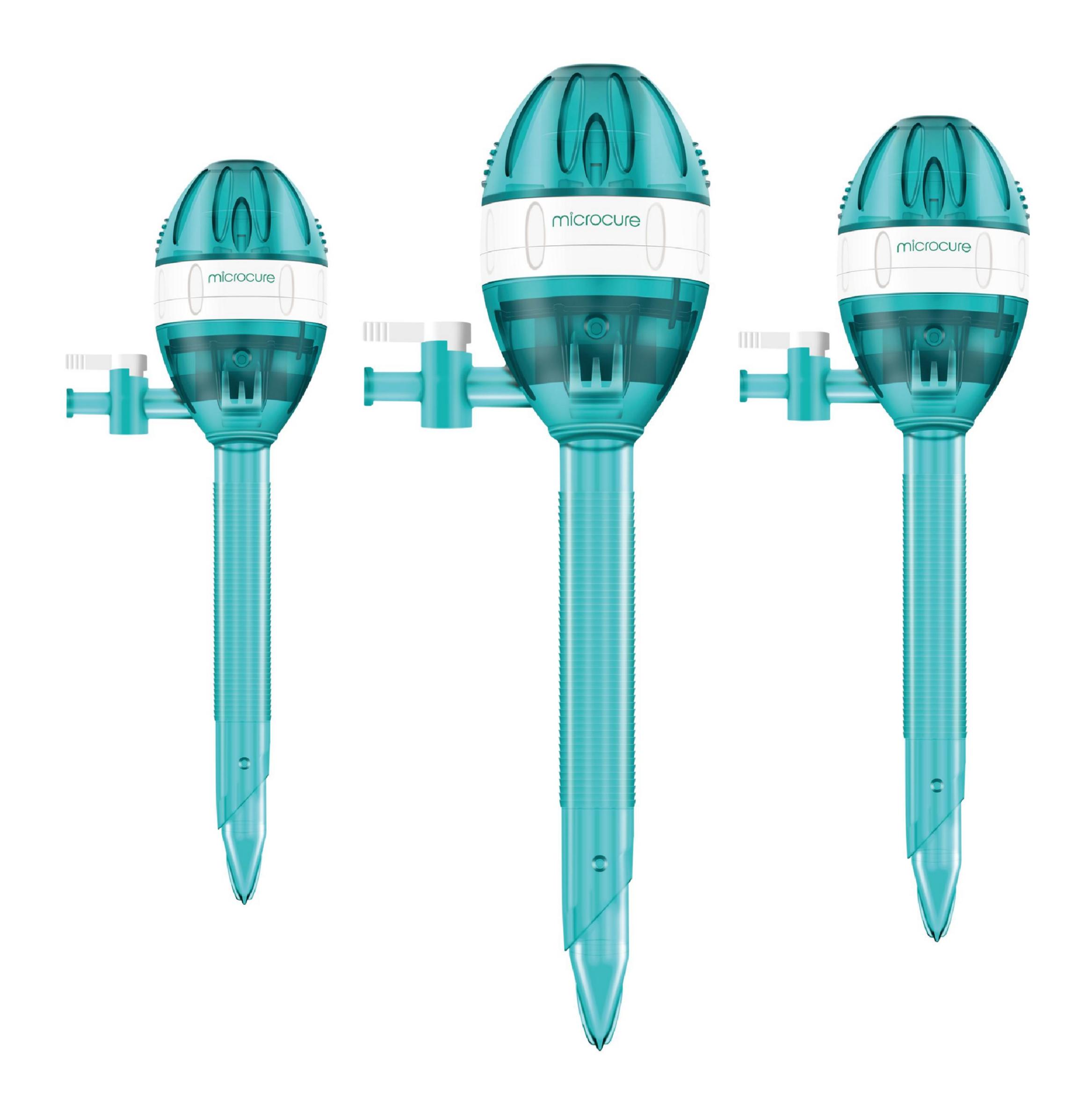

Product Description

Applicable to the establishment of microports for introduction of endoscopes and laparoscopic instruments.

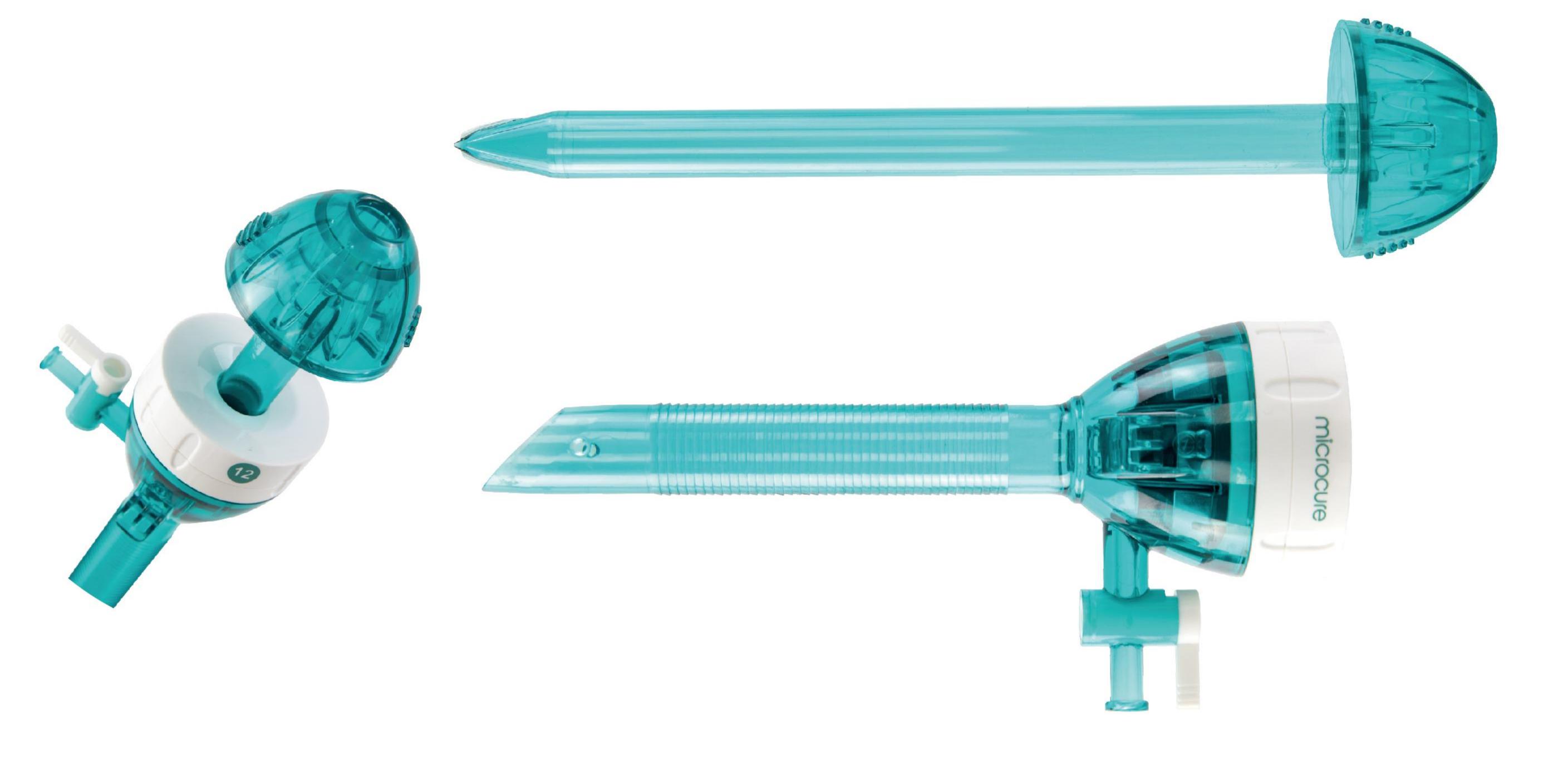

Product Advantages

01 Visible puncture

With transparent puncture rob tip,camera can be placed in puncture to prevent accidental tissue injury.


02 Bladeless obturator


Separates, instead of cuts, reducing risk of abdominal wall injuries with bleeding complications.



03 Diamond tip

Facilitate smoother insertion for easier trocar placement.

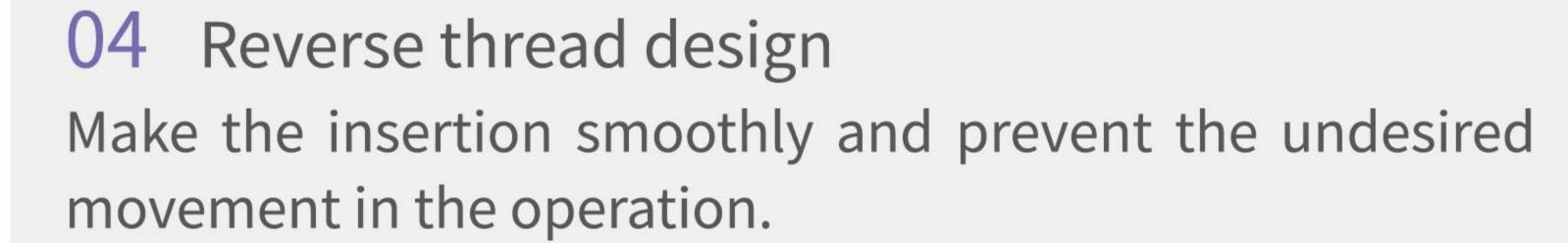

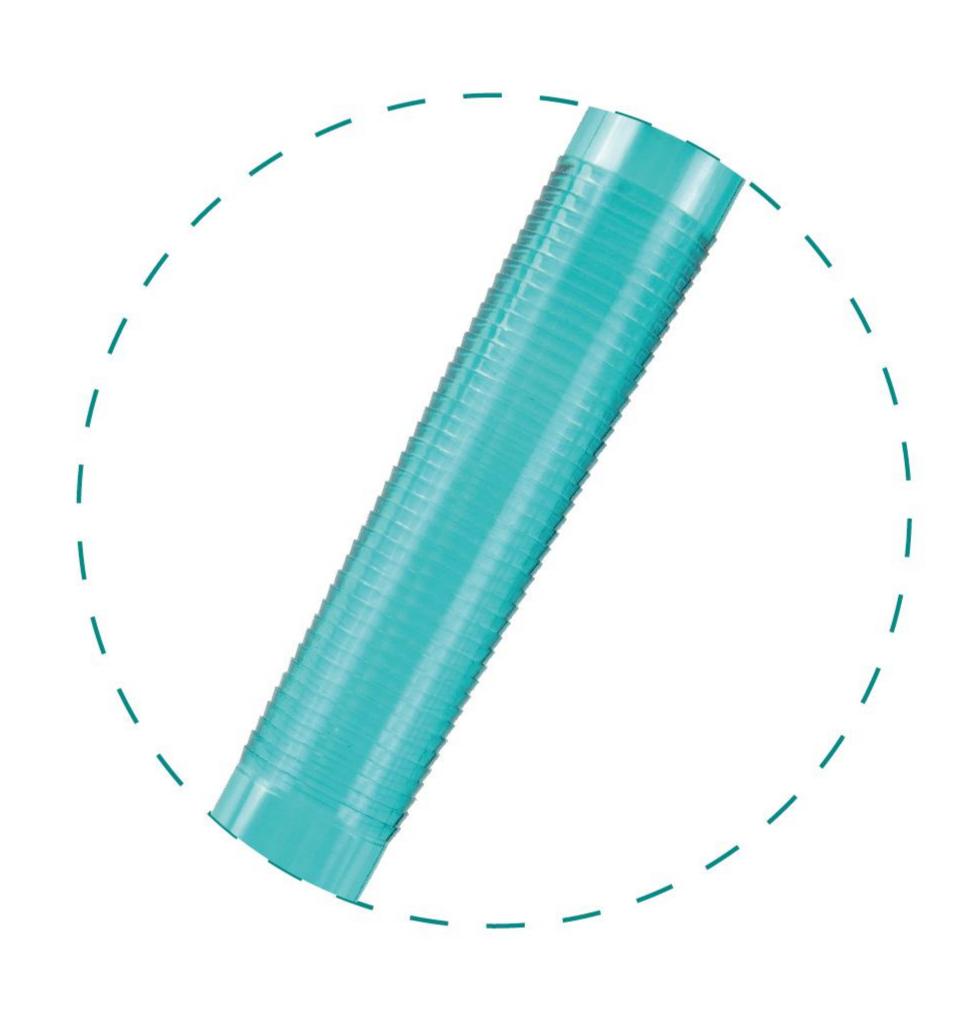



05 Removable sealing cover


The sealing cap is detachable to facilitate tissue removal during surgery.

06 Double-layer valve

The puncture cannula has better air tightness.

Double valve design



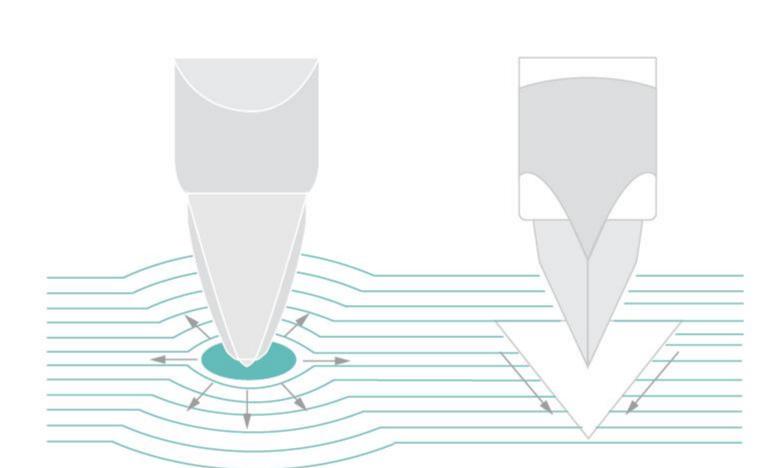


Minimally Triangular Flat knife invasive pyramid

Order Information

SPECIFICATION	MODELS	PACKAGE CONTENT	CASING WORKING LENGTH (mm)
Catanala	VT5-100	5mm Trocarx2	103
Set package	VT8-100	10mm Trocarx2	103
	VT11-100	12mm Trocarx2	103
	VT12-150	12mm Trocarx2	150
_	VT19-100	5mm Trocarx2+10mm Trocarx1	103
Set package	VT25-100	5mm Trocarx1+10mm Trocarx2	103
2	VT28-100	5mm Trocarx1+12mm Trocarx2	103
_	VT40-100	5mm Trocarx2+10mm Trocarx2	103
Set package	VT43-100	5mm Trocarx2+12mm Trocarx2	103
3	VT46-100	5mm Trocarx2+10mm Trocarx1+12mm Trocarx1	103
Set package	VT60-100	5mm Trocarx2+10mm Trocarx2+12mm Trocarx1	103
4	VT61-100	5mm Trocarx2+10mm Trocarx1+12mm Trocarx2	103
	WZDT-A-3/55	3mm Trocarx1	55
	WZDT-A-5/100	5mm Trocarx1	103
	WZDT-A-10/100	10mm Trocarx1	103
Individual	WZDT-A-12/100	12mm Trocarx1	103
	WZDT-A-15/100	15mm Trocarx1	103
	WZDT-C-5/100	5mm Trocarx1	103
	WZDT-C-10/100	10mm Trocarx1	103
	WZDT-C-12/100	12mm Trocarx1	103

*C series is bladed trocar.
*OEM services are available. microcure 124 23 microcure


DISPOSABLE TROCARS FOR THYROID SURGERY

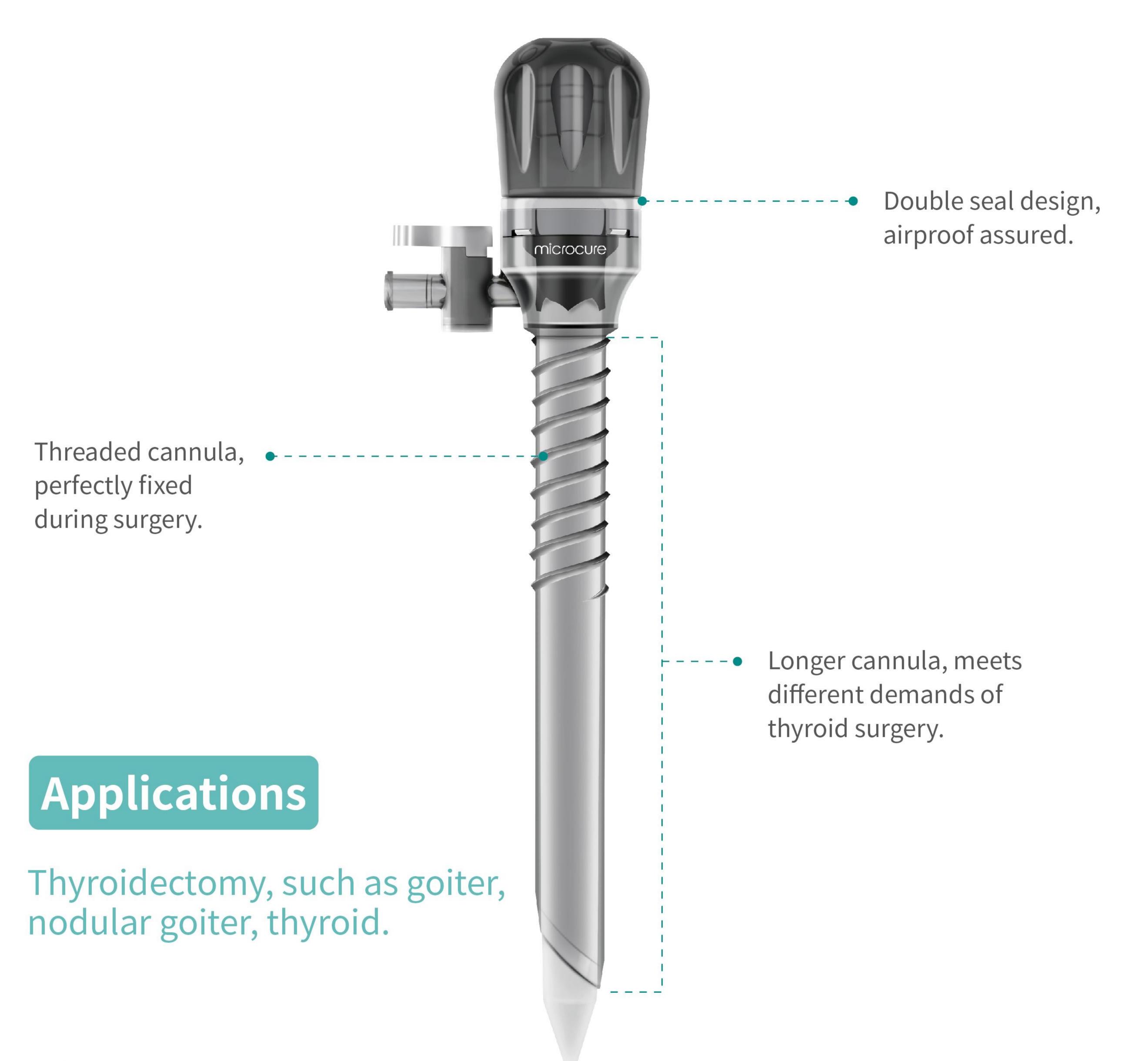

VERESS NEEDLE

Used to create a pneumoperitoneum in laparoscopic procedures


Product Advantages

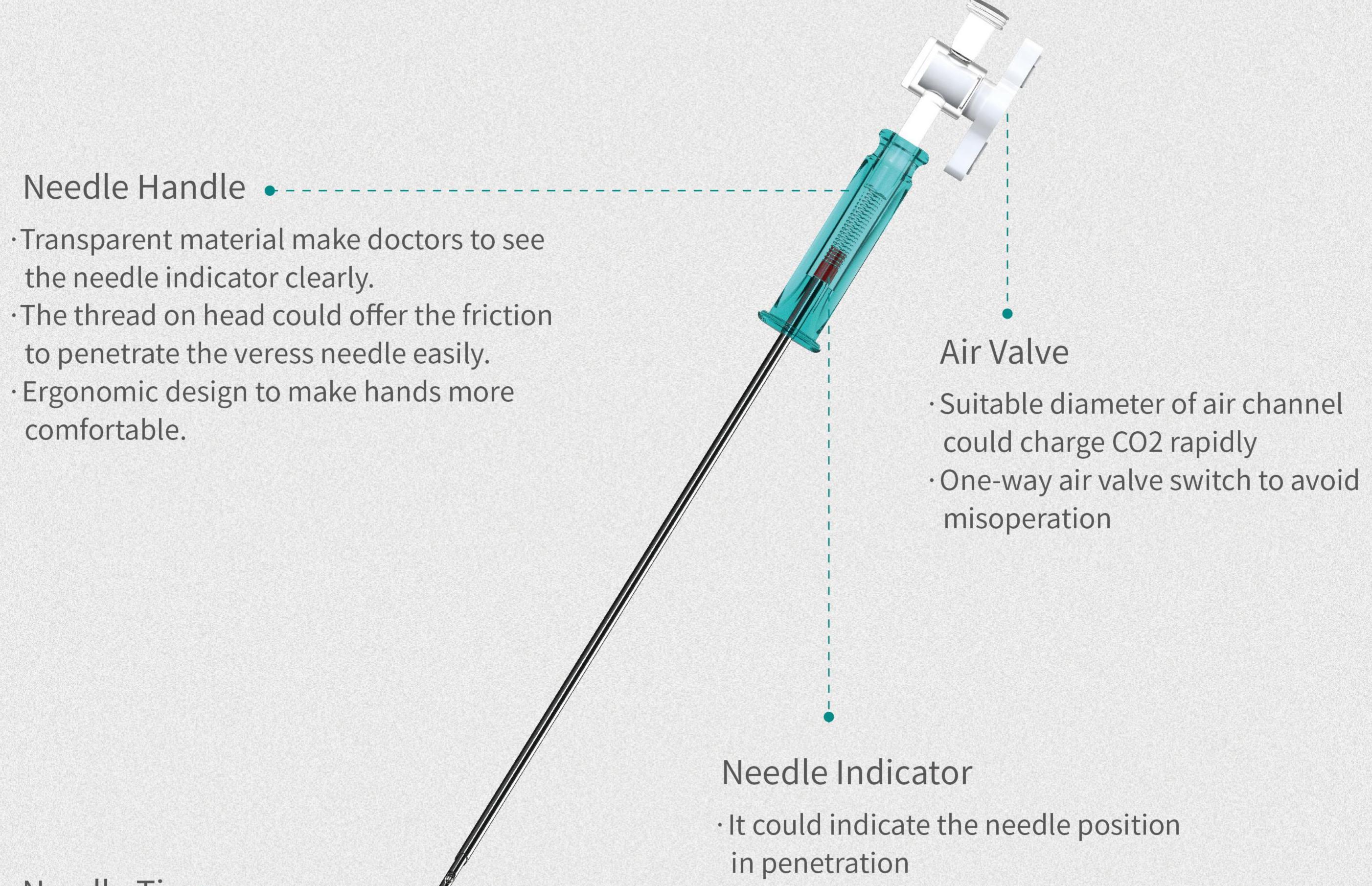

Needle Tip • -

· Blunt tip with spring loaded mecha nism to protect patients from damage. · Red color is easy to be observed

DISPOSABLE INCISION PROTECTION SLEEVE

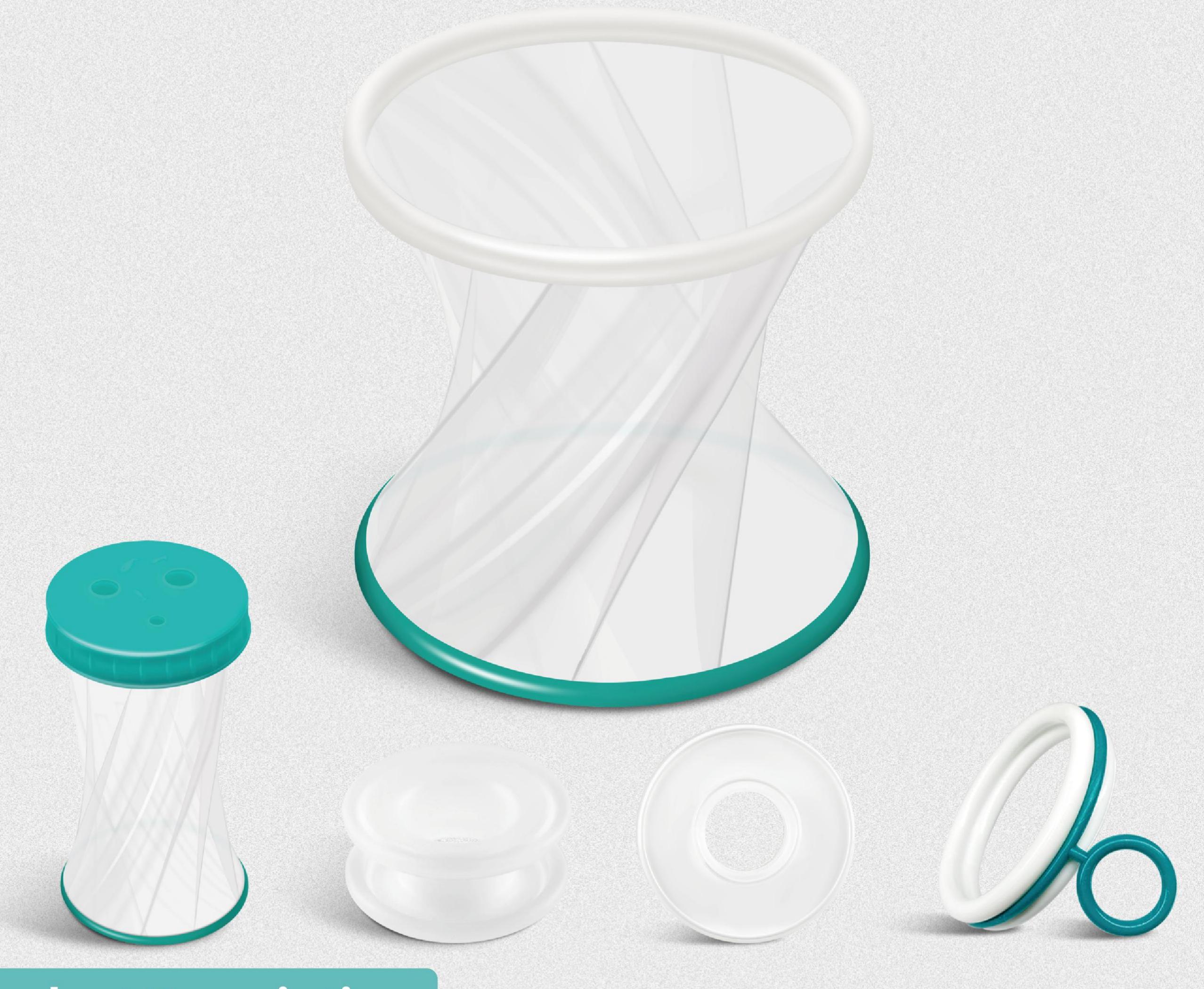

Product Description

The product series include Polymer IPS-A series and Silicone IPS-B series.

Polymer IPS-A series are made of polymer, height-adjusted. Silicone IPS-B series are made of silicone material height-fixed.

Application

Polymer IPS-A series are applicable to general surgery, G&A surgery, laparotomy and other open surgeries with large incision.

Silicone IPS-B series apply to thoracoscopic surgery, suspended laparoscopic surgery and other surgeries with small incision.

Product Advantages

01 Retraction

Provide 360° circumferential elastic retraction during procedures and evenly distribute retraction to minimize wounds and pains.

02 Protection

Protect wounds at the incision site with superior material strength, significantly avoid infection and keep incision moist.


03 Exposure

Minimize the working field for better exposure and visualization, prevent bleeding implications and free up surgeons' hands.



Product Demo


- 01 Make an appropriate sized surgical incision in the surgical area.
- 02 Insert the inner ring (Green) into abdominal cavity through the incision.
- 03 Roll the outer ring (White) inward at the position of 10 o'clock and 2 o'clock until it stops.
- 04 Through the tunnel of incision protection sleeve to operate.

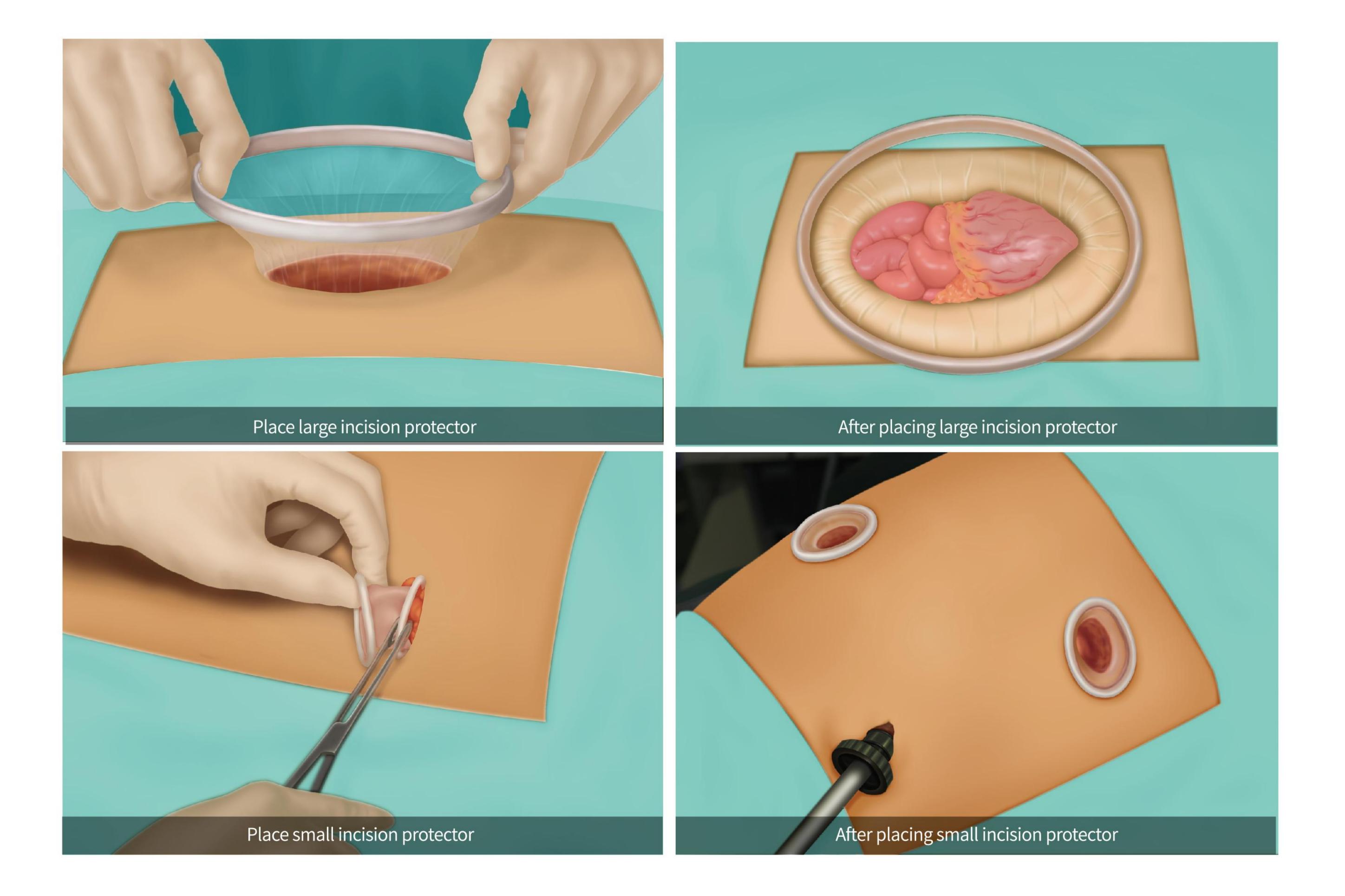

Order Information

Models		Size	Inner ring diameter (mm)	Tunnel diameter (mm)	Outer ring diameter (mm)	Tunnel length (mm)	Incision length (cm)
		50/50-50/220 (1)	50	50	50	220	1.8-3
	S	60/70-60/220 (2)	70	60	60	220	2-4
	5	60/70-60/150	70	60	60	150	2-4
		80/90-80/150	90	80	80	150	2.5-6
WZIPS-A	М	120/130-120/220	130	120	120	220	5-9
	M&L	150/160-150/220	160	150	150	220	7-12
	L	180/190-180/220	190	180	180	220	9-14
		220/230-220/250	230	220	220	250	11-17
	XL	250/260-250/100	260	250	250	100	17-20
		320/330-320/250	330	320	320	250	15-22
		40/30-15/25	40	30	15	25	1.5
WZIPS-B		50/40-25/25	40	50	25	25	2.5
	S	70/70-35/25	70	70	35	25	3.5
		70/60-45/25	60	70	45	25	1.5-2.5
		150/70-60/5	70	150	60	5	2-3

^{*(1)} and (2) are available for gastroschisis.

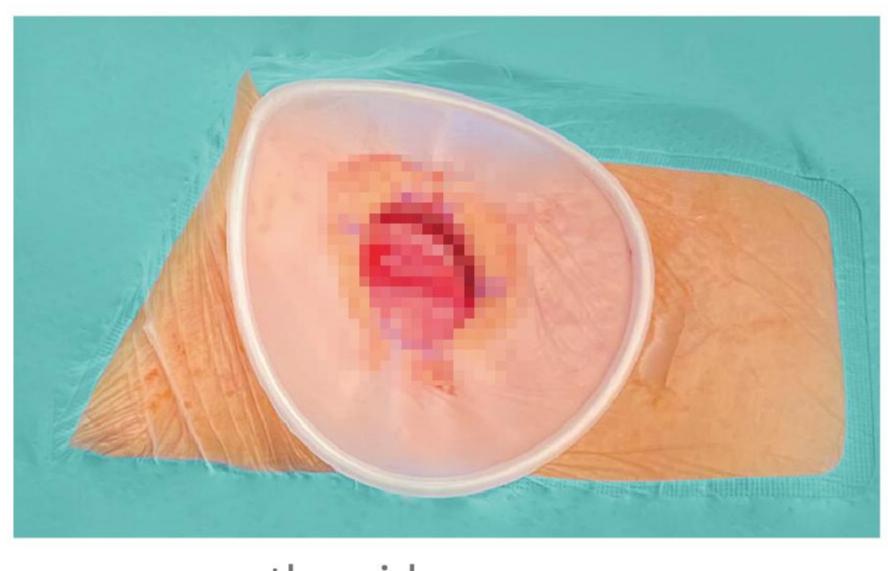

thyroid surgery



caesarean section

^{*}Our sleeves are also applicableto thyroid surgery and caesarean section.

^{*}OEM services are available.

^{*}For more model recommendations, please consult us.

SURGERY RETRACTOR SYSTEM

Product Description

The retractor system operates manually with self-locking function. It consists of a retractor ring, movable screws and retractor hooks.




Applications

Apply to temporarily separate tissues or anatomical structures to increase exposure of internal tissue or organ in surgical field.


- O1 Simple to use
 Only three steps to finish.
- 02 Meet various surgical demands
- O3 Self-retention design

 Free up surgeons' hands.
- O4 The elastic-balled hook can be more effectively fixed, prompting and not easy to slide.

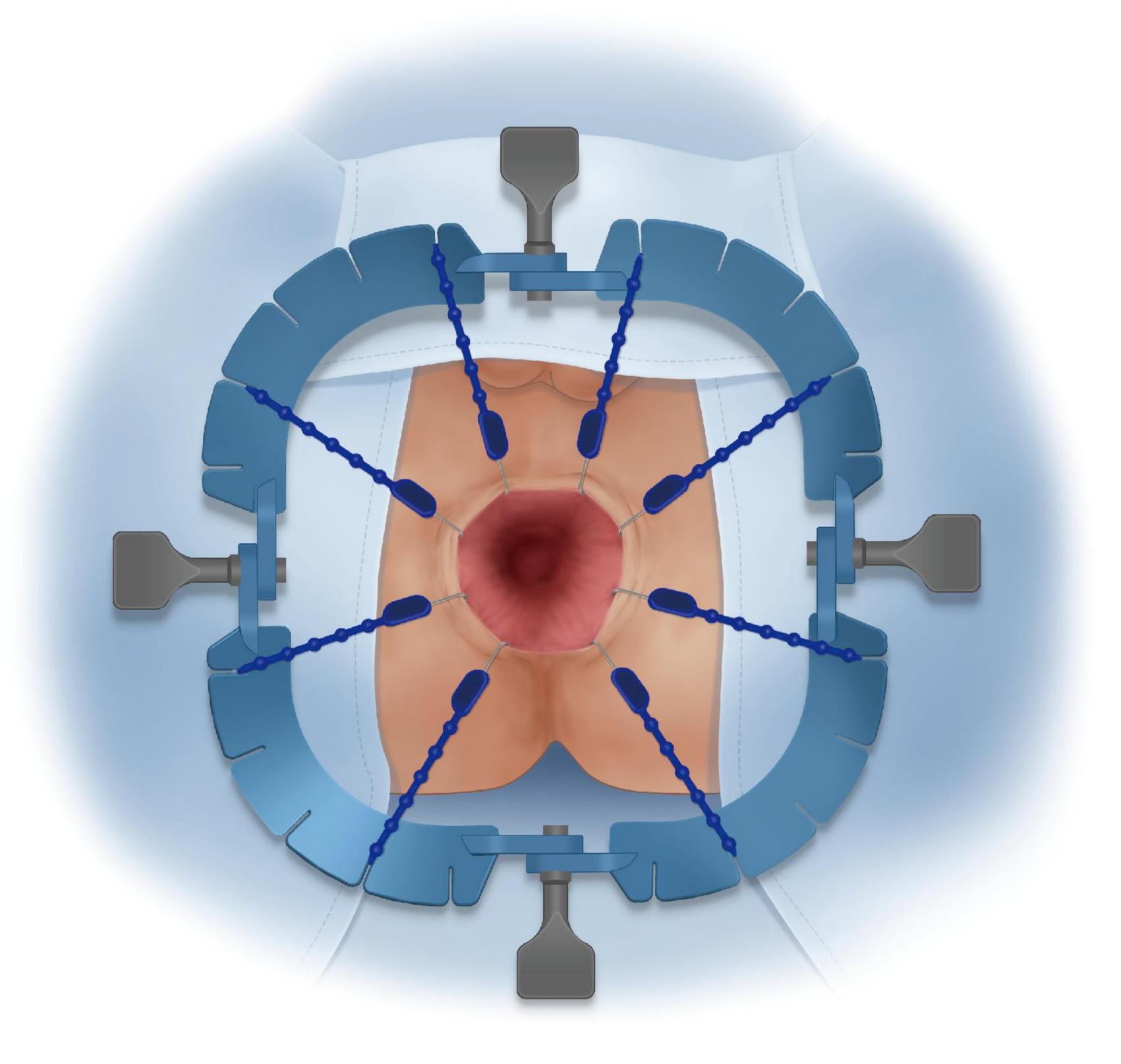



Application Departments


- 01 WZRR-A-A applicable to gynecology, head and neck surgery.
- 02 WZRR-A-B applical to otolaryngology surgery.
- 03 WZRR-A-C applical to thyroid, breast, plastic surgery.
- 04 WZRR-A-D applical to anorectal, rectal, general, gynecology.
- 05 WZRR-A-E applicable to gynecology, male urology.
- 06 WARR-A-F applical to thyroid, urology surgery.
- 07 WZRR-A-G applicable to general surgery, head and neck surgery.
- 08 WZRR-A-H applical to gynecology, urology surgery.

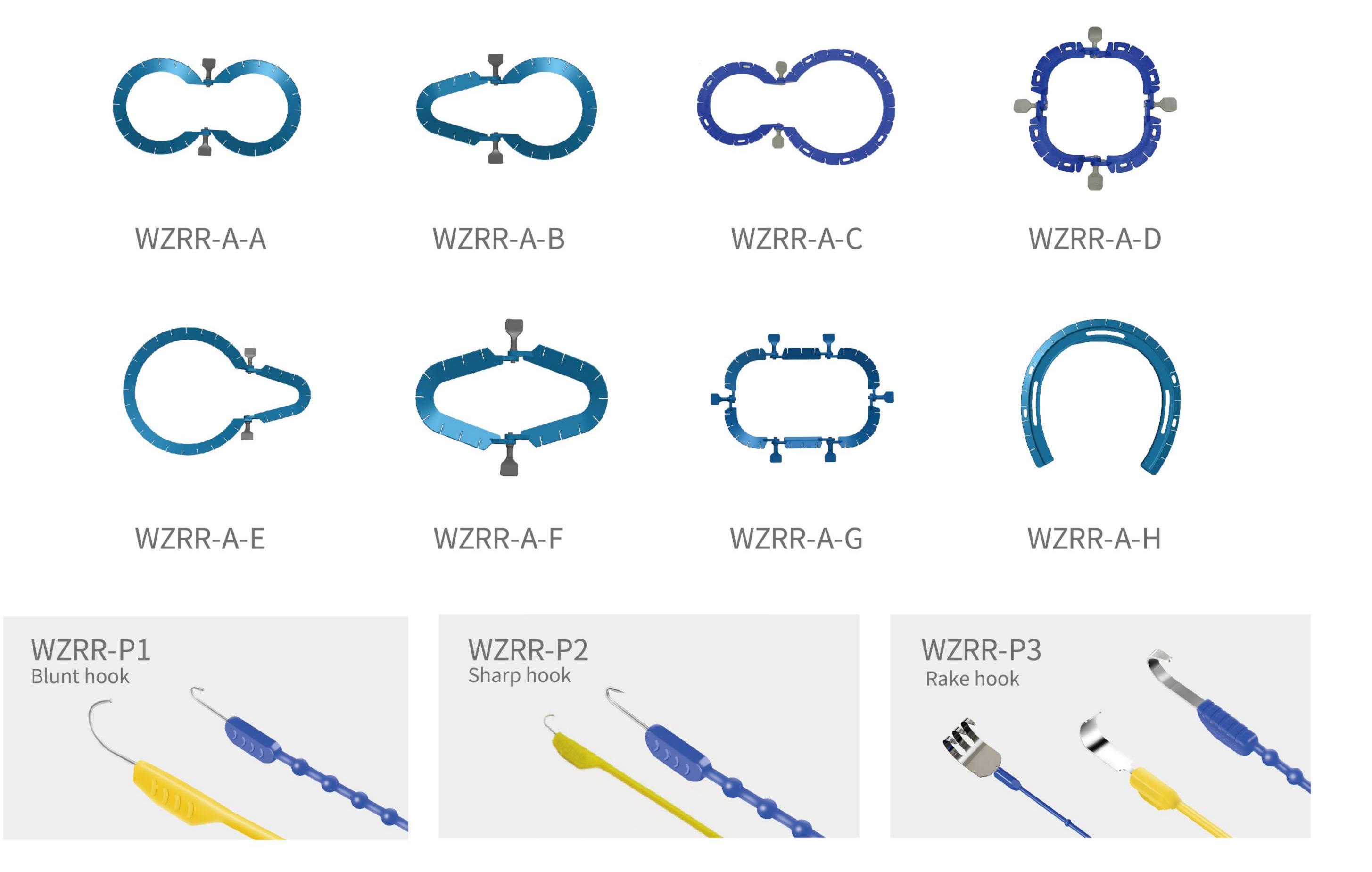

Product Demo

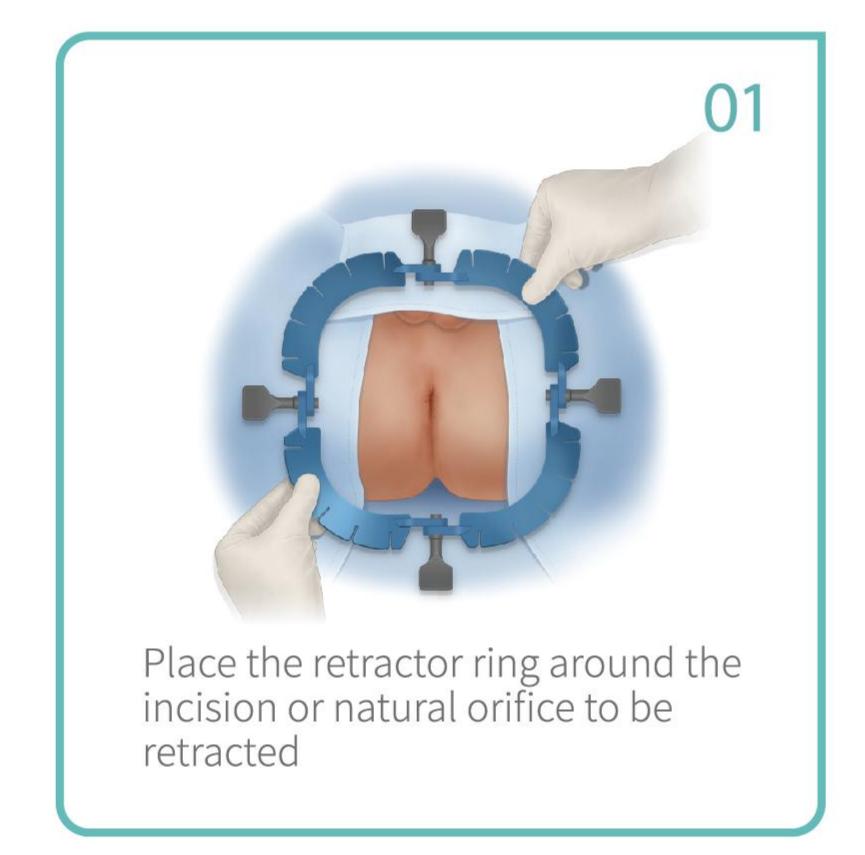

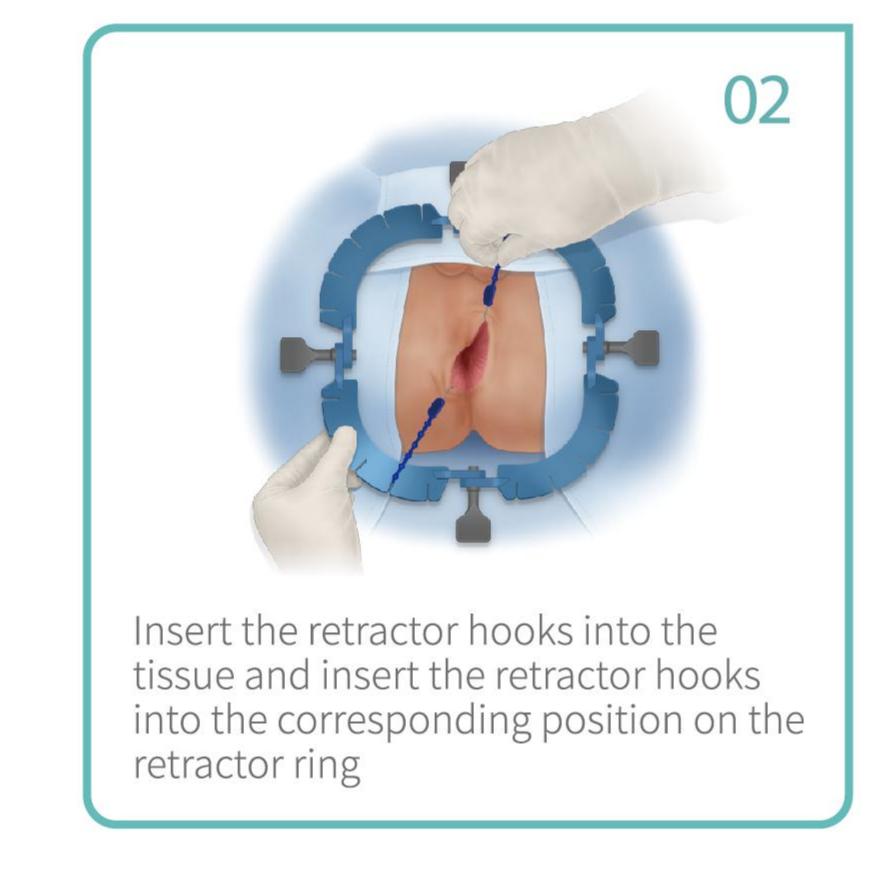

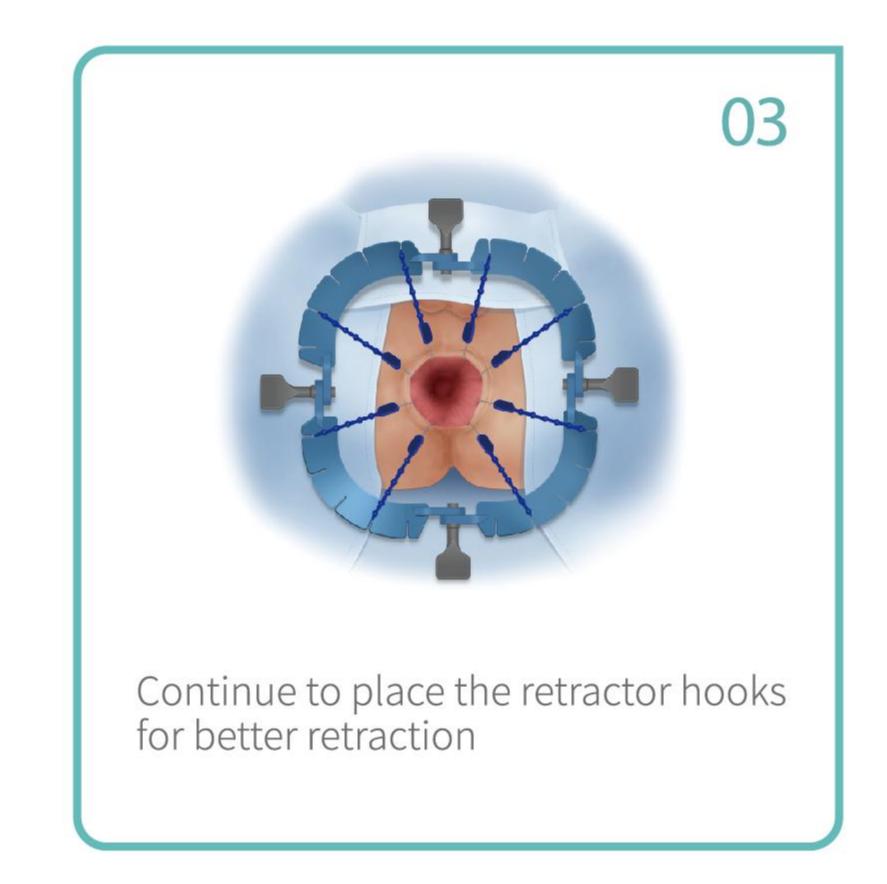

Order Information

PRODUCTS	MODELS
Retractor Ring	WZRR-A-A
	WZRR-A-B
	WZRR-A-C
	WZRR-A-D
	WZRR-A-E
	WZRR-A-F
	WZRR-A-G
	WZRR-A-H
Retractor Hook	WZRR-P1-S-1
	WZRR-P1-S-2
	WZRR-P1-L-2
	WZRR-P2-S-1
	WZRR-P3-L-1
	WZRR-P1-L-1
	WZRR-P2-L-1
	WZRR-P3-L-3
	WZRR-P3-L-4

*In the model "WZRR-P1-L-1", "P1" indicates hook's type, "L-1" indicates retractor's size and the number of claws.

*OEM services are available.

33 | microcure | 34

DISPOSABLE SAMPLE RETRIEVAL BAG

Product Advantages



01 Detachable handle

Convenient to separate the handle and rob and ergonomically designed handle makes it comfortable for pushing and pulling.



02 Material

The TPU material bag is soft and transparent. Ni-Ti metallic rim enables bag to open automatically,quickly and completely.



03 Retrieval bag

The retrieval bag can be opened automatically, quickly and completely and it also can be opened closed repeatedly.



Applications

Procedural Applications Facilitates a wide range of laparoscopic procedures: Cholecystectomy, Oophorectomy, Appendectomy, Nephrectomy, Splenectomy, Sleeve, Gastrectomy.

APPLICABLE DEPARTMENT

GENERAL SURGERY

UROLOGY SURGERY

GYNAECOLOGY SURGERY

THORACIC SURGERY

Order Information

MODELS	CAPACITY(ml)	DIAMETER(mm)	SHORT AXIS(mm)	LONG AXIS(mm)
WZSRB-A-40	35	40		N/A
WZSRB-A-60	80	60	N/A	
WZSRB-B-80	200		82	92
WZSRB-B-100	350	N/A	103	115
WZSRB-B-130	600		118	154


^{*}disposable sample retrieval bag is divided into a circular bag opening (A) and an oval bag opening (B).
*We also provide re-openable 350ml retrieval bag.

^{*}OEM services are available.


DISPOSABLE LIGATION LOOP


Applications


It is suitable for ligation of lumen-like tissues (such as appendix, uterine fibroids, etc.) in various endoscopic operations.

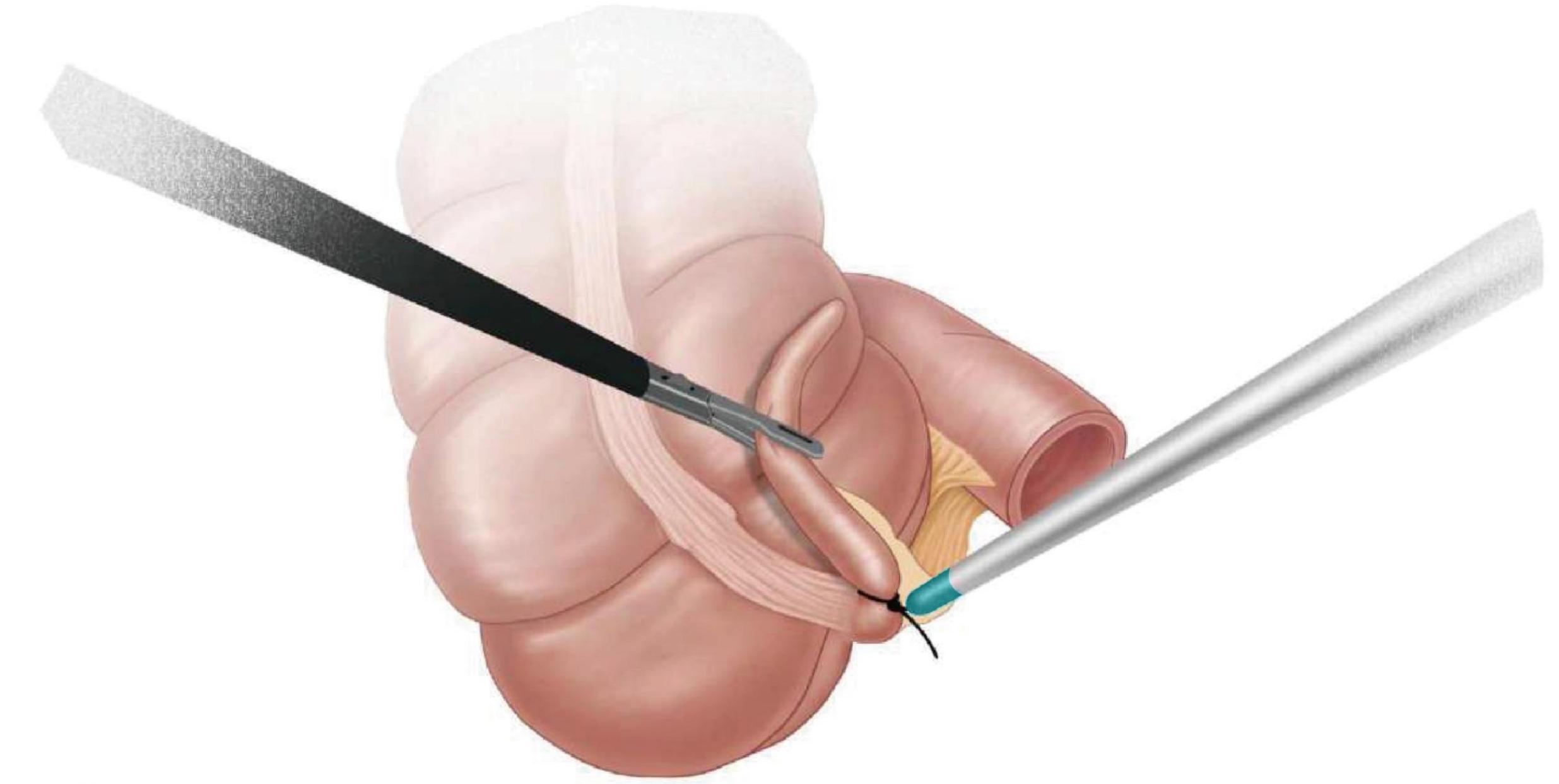


Product Demo

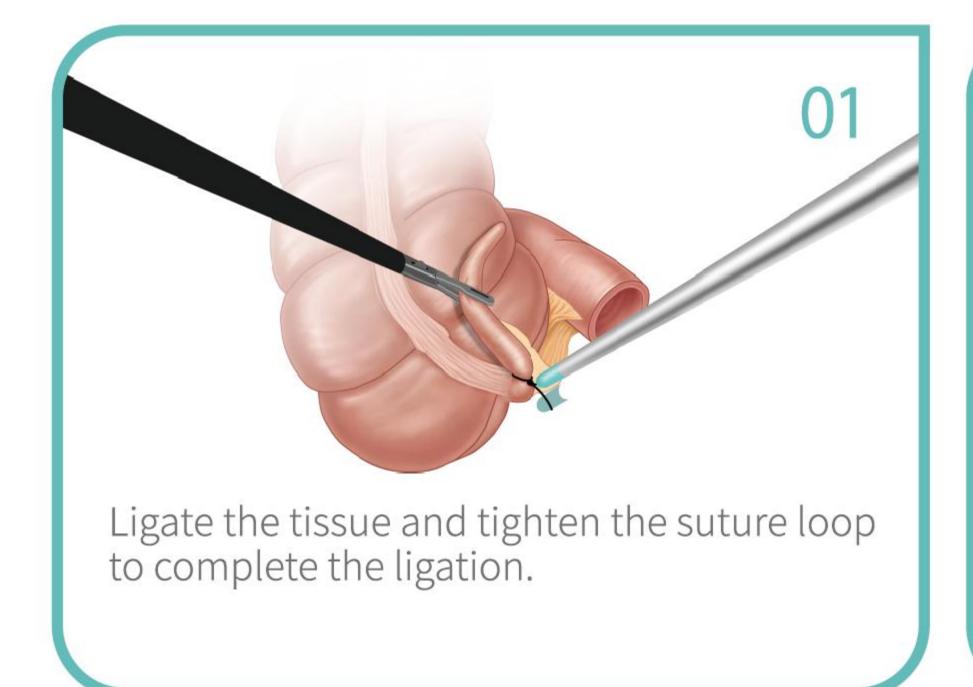

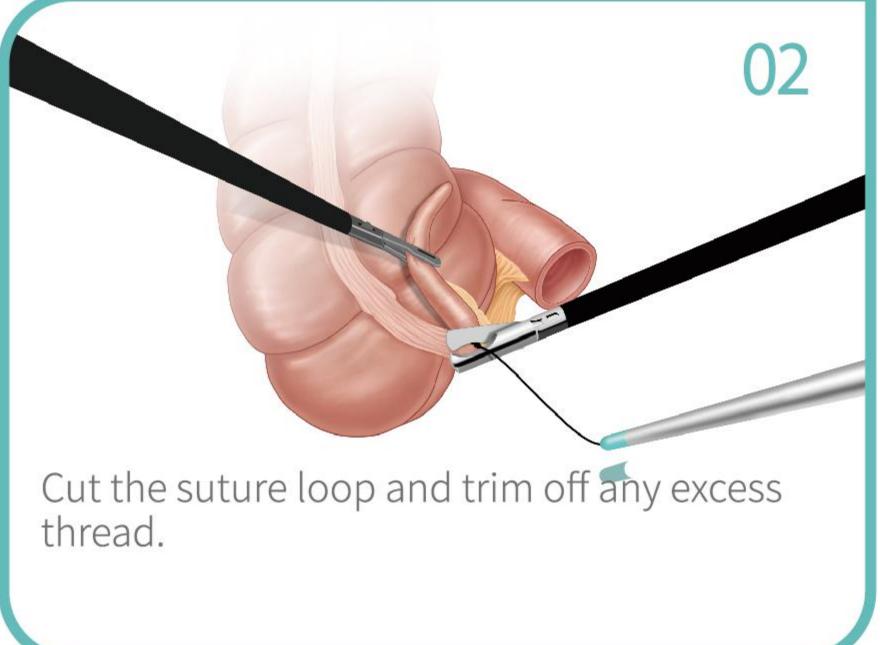

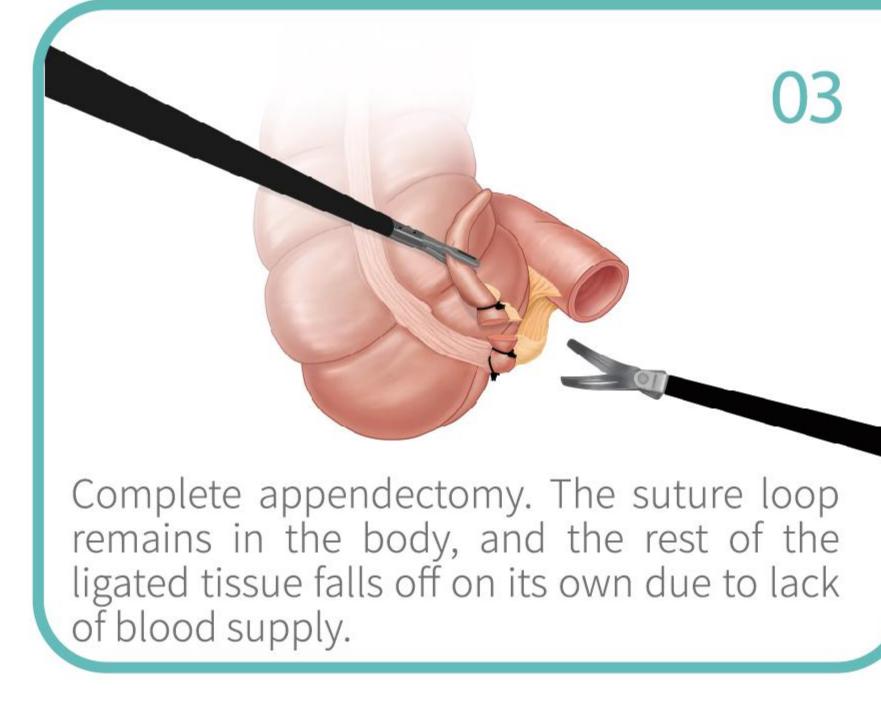

Product Advantages

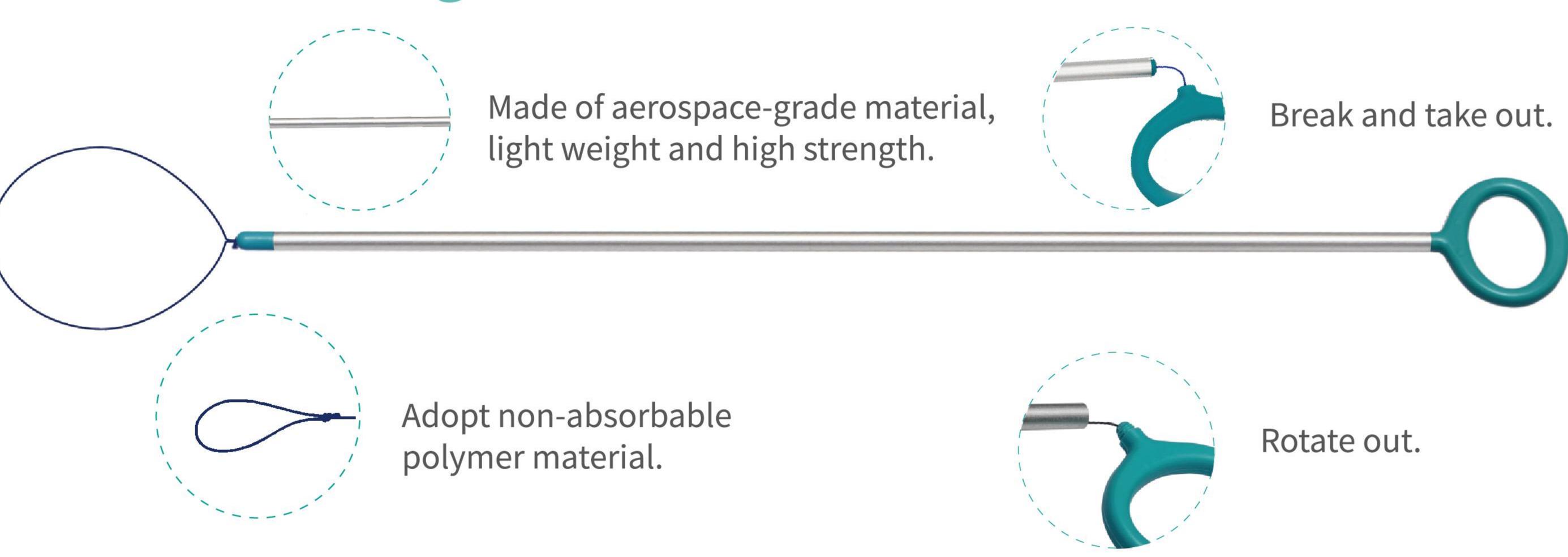

01 Handle with two types of usage

unscrew or break off directly.

02 Connecting rod

made of aluminium alloy with high strength and lower weight.

03 Ligature


the non-absorbale surgical suture has steady physicochemical proporties which decreases the risk of infection.


Order Information

MODELS	TOTAL LENGTH(mm)	LIGATURE MODELS	LIGATURE DIAMETERS(mm)
WZDLL-A-345/1	345	1	0.400-0.499
WZDLL-A-345/0	345	0	0.350-0.399
WZDLL-A-345/2-0	345	2-0	0.300-0.349

*OEM services are available.


Certificate

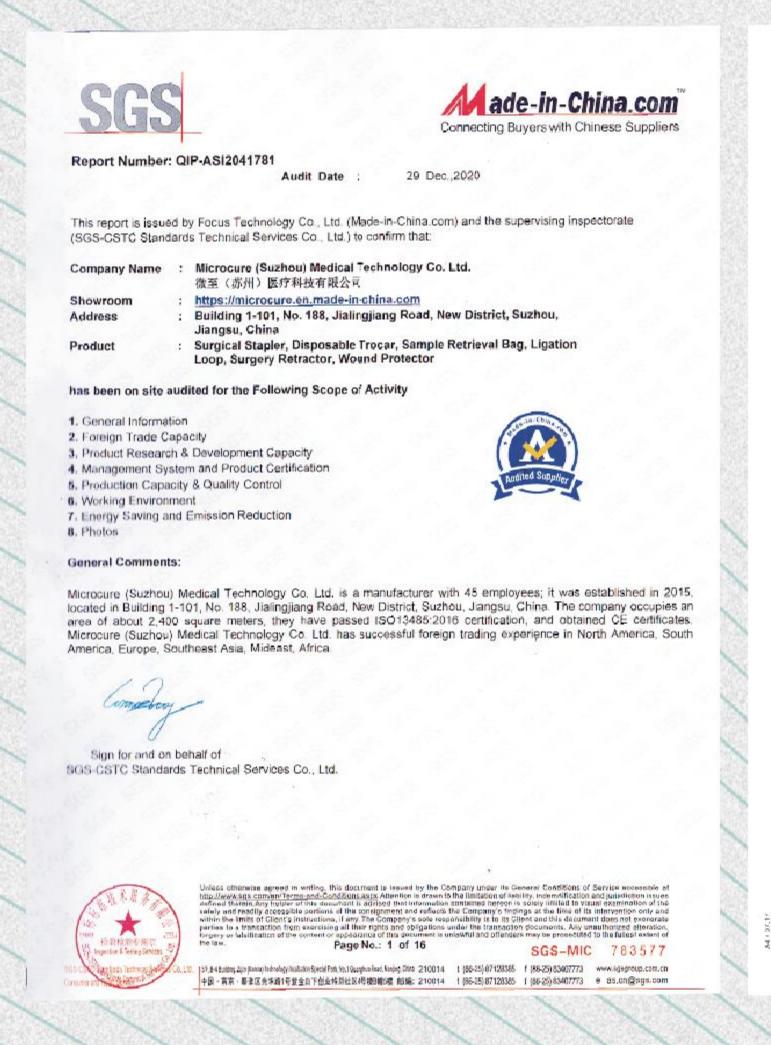

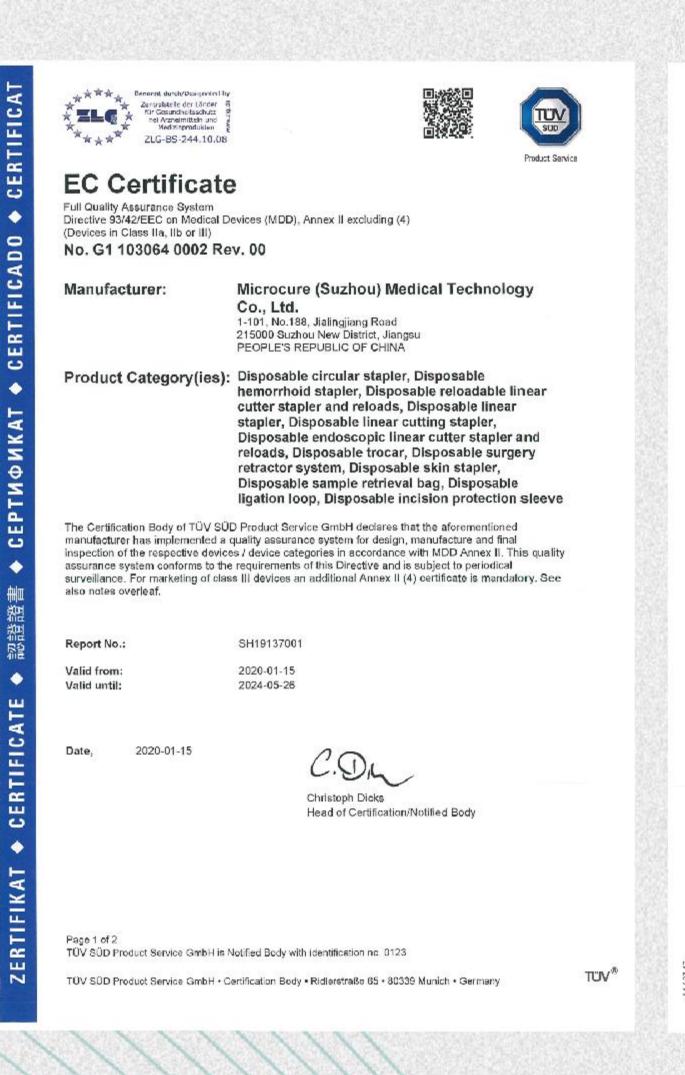



In perfection we believe

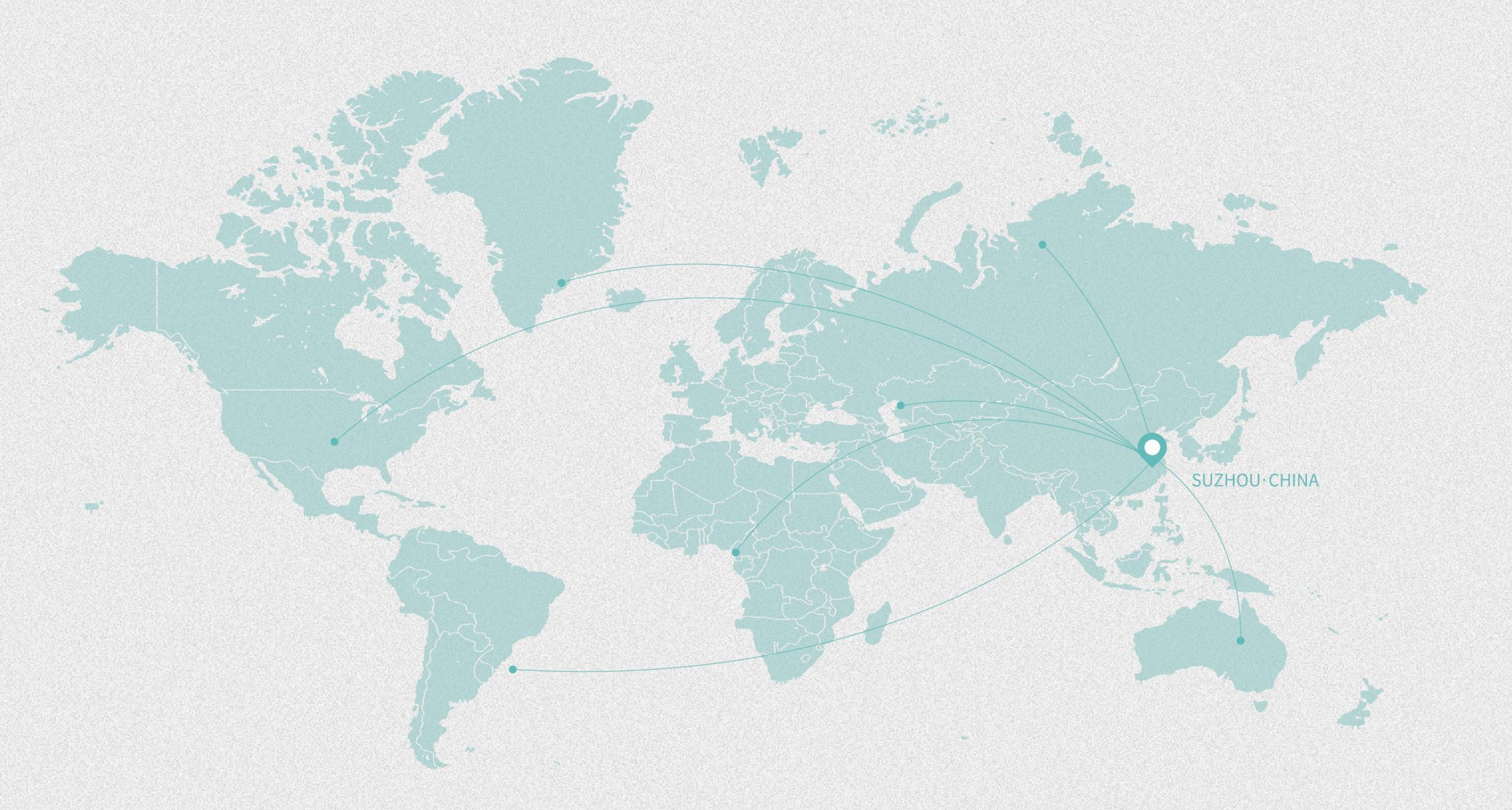

Michael