



# One system. Many options.



#### Standard or over-the-wire delivery

The only conical retrievable IVC filter cleared by the FDA for over-the-wire delivery.

- The unique filter design allows for delivery over a wire ≤ .035" using the standard kit components
- The wire acts as a guide to deploy the filter straight and in the center of the IVC positioning it for easier retrieval
- Over-the-wire delivery provides stability to the apex during deployment for precise and predictable positioning in IVC

#### More access for delivery


The only IVC filter cleared by the FDA for delivery via jugular, femoral, antecubital, and popliteal access.

 Save time while providing PE protection during DVT treatment by delivering the OptionELITE via popliteal access to avoid repositioning and redraping the patient

# For temporary or permanent placement

- Unique anchor and apex design makes retrieval easy when the need for PE protection subsides
- A robust design with an established positive safety profile provides confidence for long-term PE protection

## Option ELITE

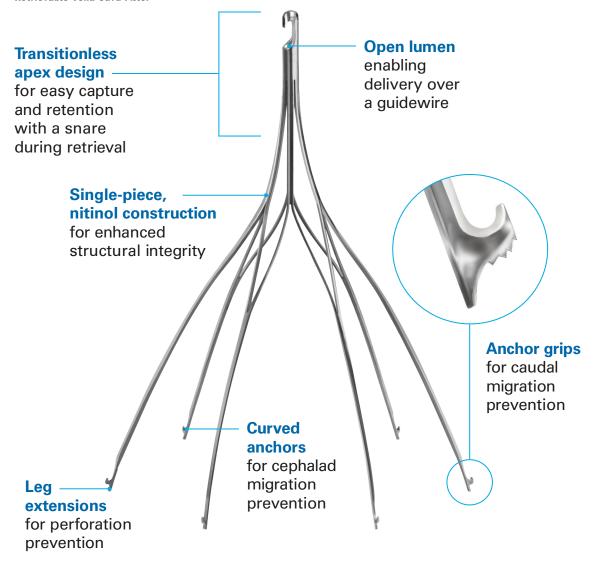

Initial Clinical Trial Results <sup>1</sup> (100 Patients)
2.0% Migration
0.0% Fractures
1.0% Perforations
3.0% Occlusion

IVC Filter Performance Data <sup>2</sup> (165 Patients)
0.0% Migration
0.0% Fractures
1.2% Perforations
0.6% Occlusion

With over a decade of positive, consistent, and proven clinical performance, the Option design has achieved 100% technical and retrieval success with O fractures and O migrations in a 165-patient study.<sup>2</sup>

### Components

### **Curved Pusher Delivery Sheath Dilator Hemostatic** valve to prevent **Deployment mark** bleed back for alignment during delivery Integrated with an over the side port wire technique for hand injection **Deployment mark** for alignment during delivery with a standard technique Low-profile design with a 5 FR ID

Curved pusher aids centering and resists tilt during standard delivery



OptionELITE retrievable vena cava filter



Universal color-coded filter cartridge for placement via four different access sites



Valve cap maintains hemostasis during post-placement cavagrams

Radiopaque makers spaced 32mm apart

### **Ordering Information**

OptionELITE Inferior Vena Cava Filter System						
CATALOG #	DESCRIPTION	CAVA DIAMETER	SHEATH LENGTH	SHEATH SIZE		
352506070E	Femoral or Jugular Access	up to 30 mm	70 cm	5 F		
352506100E	Antecubital or Popliteal Access	up to 30 mm	100 cm	5 F		

Single-Loop Retrieval Kit						
CATALOG #	INNER SHEATH	OUTER SHEATH	SNARE CATHETER			
380911120	9F (ID) x 64 cm	11F (ID) x 60 cm	7F (OD) x 76 cm			

Triple-Loop Retrieval Kit						
CATALOG #	INNER SHEATH	OUTER SHEATH	SNARE CATHETER			
380911330	9F (ID) x 64 cm	11F (ID) x 60 cm	7F (OD) x 76 cm			

V·Stick Vascular Access Sets with Co-axial Introducer (Box of 10)							
	SHEATH	NEEDLE		DILATOR		GUIDEWIRE	
CATALOG #	SIZE	GUAGE x LENGTH	ECHOGENIC	STIFFNESS	NUT COLOR	DIAMETER x LENGTH	MATERIAL
390204400E	4 F	21G x 3.8 cm	No	Standard	Red	.018" x 40 cm	Stainless Steel
391504000E	4 F	21G x 7 cm	Yes	Standard	Red	.018" x 40 cm	Stainless Steel
391504300E	4 F	21G x 7 cm	Yes	Standard	Red	.018" x 40 cm	Nitinol
391504400E	4 F	21G x 7 cm	No	Standard	Red	.018" x 40 cm	Nitinol
391504500E	4 F	21G x 7 cm	No	Standard	Red	.018" x 40 cm	Stainless Steel
391505000E	5 F	21G x 7 cm	Yes	Standard	Red	.018" x 40 cm	Stainless Steel
391505300E	5 F	21G x 7 cm	Yes	Standard	Red	.018" x 40 cm	Nitinol
391505400E	5 F	21G x 7 cm	No	Standard	Red	.018" x 40 cm	Nitinol
391505500E	5 F	21G x 7 cm	No	Standard	Red	.018" x 40 cm	Stainless Steel
391604000E	4 F	21G x 7 cm	Yes	Stiff	White	.018" x 40 cm	Stainless Steel
391604300E	4 F	21G x 7 cm	Yes	Stiff	White	.018" x 40 cm	Nitinol
391604400E	4 F	21G x 7 cm	No	Stiff	White	.018" x 40 cm	Nitinol
391604500E	4 F	21G x 7 cm	No	Stiff	White	.018" x 40 cm	Stainless Steel
391605000E	5 F	21G x 7 cm	Yes	Stiff	White	.018" x 40 cm	Stainless Steel
391605300E	5 F	21G x 7 cm	Yes	Stiff	White	.018" x 40 cm	Nitinol
391605400E	5 F	21G x 7 cm	No	Stiff	White	.018" x 40 cm	Nitinol
391605500E	5 F	21G x 7 cm	No	Stiff	White	.018" x 40 cm	Stainless Steel

Guidewires (Box of 10)							
CATALOG #	DESCRIPTION	DIAMETER	LENGTH	TIP	MATERIAL		
395369	Fixed Core with PTFE coating	035 in	260 cm	Straight	Stainless Steel		
395370	Fixed Core with PTFE coating	038 in	260 cm	Straight	Stainless Steel		
395722	Fixed Core with PTFE coating	035 in	200 cm	3mm J-Tip	Stainless Steel		
395794	Fixed Core with PTFE coating	038 in	260 cm	3mm J-Tip	Stainless Steel		
395170	Bentson Wire with PTFE coating	035 in	180 cm	Straight	Stainless Steel		
390033	Bentson Wire with PTFE and Heparin coating	035 in	260 cm	Straight	380911120		

<sup>© 2020</sup> Argon Medical Devices, Inc. All rights reserved. Argon Medical Devices, the stylized logo, and the product and service marks mentioned herein are trademarks or registered trademarks of Argon Medical Devices, Inc. in the United States and other countries. All other trademarks are the property of their respective owners.



<sup>1.</sup> The Safety and Effectiveness of the Retrievable Option™ Vena Cava Filter: A United States, Prospective Multicenter Clinical Study Matthew S. Johnson, MD, James F. Benenati, MD, John A. Kaufman, MD, et al. J Vasc Interv Radiol 2010; 21:1173–1184.

<sup>2.</sup> Iqbal SI, Elmi A, Hedgire S, et al. Preliminary Experience With Option Inferior Vena Cava Filter. Vasc and Endo Surg 2012; 47(I) 24-29.