

# FemoSeal™

## Vascular Closure System




FemoSeal is a medical device indicated for use in closing and reducing time to hemostasis at the femoral arterial puncture site in patients who have undergone diagnostic or interventional procedures with 7F (2,33 mm) procedural sheath or smaller (minimal artery diameter of 5 mm).

### Product Characteristics

- FemoSeal consists of two fully bioabsorbable polymer discs held together by bioabsorbable sutures. Hemostasis is achieved by mechanically sealing the arteriotomy between the two discs
- Single use, one year shelf life after sterilization
- Sterilized by ethylene oxide. Do not re-sterilize
- Available in select markets
- 10 units per box

#### Storage, packaging and disposal:

- FemoSeal should be stored in a cool location (room temperature not exceeding 25 °C)
- FemoSeal contains resorbable materials that degrade by exposure to heat and moisture; therefore, the device may not be resterilized
- FemoSeal is sterile and non-pyrogenic in unopened and undamaged package
- Dispose of contaminated units, components, and packaging materials utilizing standard hospital procedures and universal precautions for biohazardous waste

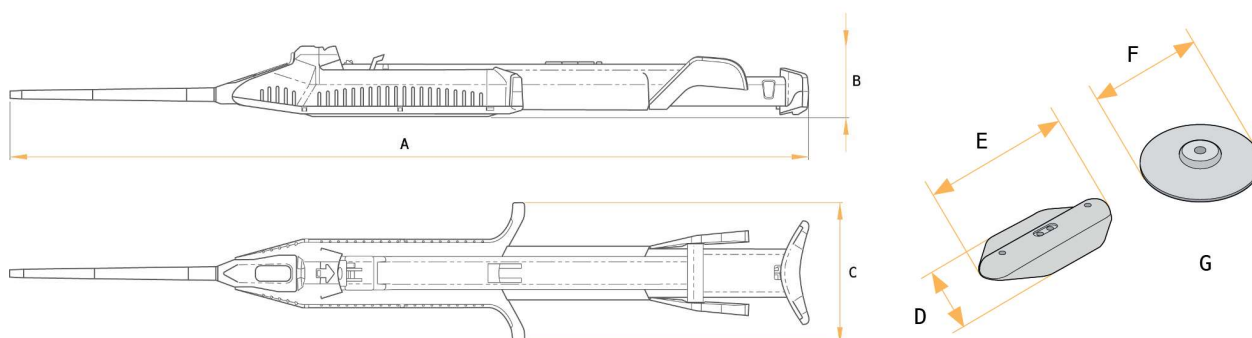

A 282 mm  
B 25 mm  
C 50 mm  
D 5 mm  
E 10 mm

F 5 mm  
G 1 mm thick

### General specifications

Effective Length Sheath	70 mm
French Size	7 Fr
Guidewire - Diameter	0.038 in / 0.97 mm

### Ordering information

C11202

Please quote above item reference codes when placing an order

Unit per box: 10 pcs.

Please quote above item reference codes when placing an order