
JJJ DAHLHAUSEN[®]

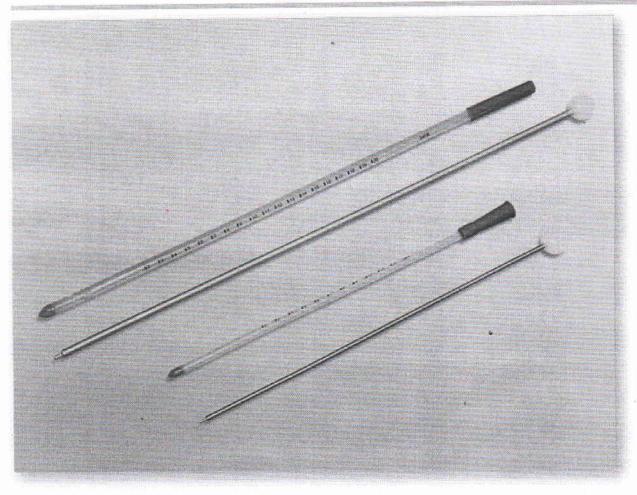

Trocar Catheter

The Trocar Catheter is introduced into the pleural cavity after puncture incision. The light metal mandrin is removed and the remaining catheter is connected to the thorax drainage unit.

The Trokar catheter is graduated in steps of 1 cm starting beneath the last lateral eye, for a safe evaluation of the position. The Trocar catheter has a closed atraumatic tip which reduces injuries.

Sizes 10, 12 and 16 are equipped with a funnel-shaped connector.

X-ray contrastable:	X-ray contrast thread	
Scale value:	1 cm	
Sterile:	Yes	

Art.no.	Outer diameter	Length	Number (eyes)
0200030010	⊜ 10 Fg	24 cm	2
0200030012	12 Fg	24 cm	2
0200030016	16 Fg	26 cm	3 .
0200030020	20 Fg	40 cm	3
0200030024	24 Fg	40 cm	3
0200030028	28 Fg	40 cm	3
0200030032	32 Fg	40 cm	3