

Specifications

Power Supply	100~240VAC 50/60Hz 200VA
Flow Rate	0.1~10ml/s
Pressure Limit	325PSI
Auto Filling Rate	8ml/s
Air Purge Rate	1.0~9.9ml/s
Piston Rod Speed	9.9ml/s
Phases	6
Injection Programs	2000
Injection Records	2000
Volume	A: 1~200ml B: 1~200ml
Multiple Language User Interface	

High Quality Disposables



200ml Syringe for Contrast
200ml Syringe for Saline

antmed

Shenzhen Antmed Co., Ltd.

18 Jinhui Ave., Pingshan New District, Shenzhen 518122, China

Tel: +86 755 8606 0992 Fax: +86 755 8606 0991

E-mail: info@antmed.com <http://www.antmed.com>



CE 0123

P/N: EN-CT(D-L)-202005

innaStar



CT Dual Injector

Leading Contrast Media Injectors Forward

antmed
Partner For Life

CT Dual Injector

Introducing an all-new ImaStar CT dual injector, offering an extended list of product features, reliable performance, and a safe patient experience.



Ease of Operation and Clinical Efficiency

- Large touch-screen display on the injector head provides convenient access for loading, filling and purging air during set up
- Optional IV pole
- Test injection mode

Independent Dual Motor Drive System

- Dual motor drive capability for time saving loading and filling of syringes
- Allows independent injection of contrast and saline

Dual Head Flexibility

- Simultaneous injection protocol designed for advanced cardiac CTA and Spectral CT imaging exams
- Optimization of contrast usage
- Post injection saline flush separate from contrast injection to reduce contrast wastage

Simultaneous Flow Feature

- Allow simultaneously filling, air purging, contrast and saline injection
- Allows the control of contrast density in both sides of the ventricles
- Significantly reduces CT artifact
- Able to observe the right coronary arteries, right atrium and right ventricle in the same image

Heater

- Ensures contrast maintains at ideal temperature

Color Touch-Screen

- 12.1 and 7" color touch-screens for easy operation
- Clear and concise information display to ensure user operational confidence

Injection Record Storage

- Stores up to 2000 injection history records

Auto Detect

- Syringe detection feature automates air purging and piston rod operation

Connect to Other Imaging Equipment



- Connect with other imaging equipment for enhanced imaging, use contrast more efficiently and minimize radiation
- Option

A Variety of Configuration Options

- Stores up to 2000 customized injection programs
- Program storage capacity can be increased based on user needs

Real-time Pressure Curve

- Display provide real-time pressure monitoring curve

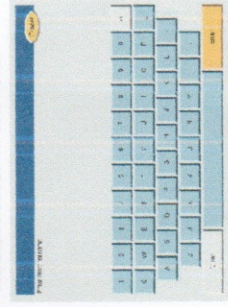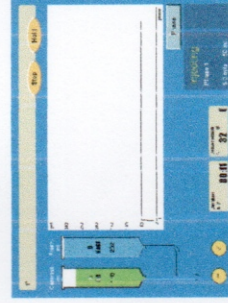

Safety

- To ensure the safety of the injection, there is a double confirmation feature with separate prompts to confirm air purging and injection ready status before any injection
- Automatically slows down the injection when approaching to pressure limit
- Injection stops when pressure limit is reached. Alerts with both light and sound
- Waterproof injector head designed to minimize damage from contrast and saline leakage
- Ability to pause and resume while injection is underway

Bluetooth communication