

Alterna[®] Free

*Spiral security with
more freedom*

Alterna[®] Free


Alterna[®] Spiral adhesive

The spiral adhesive is unique to Alterna[®].
The spiral form helps to ensure that
Alterna[®] is both secure and skin-friendly.

NEW
Increased
flexibility

Skin-friendly spiral, keeps
the skin comfortable and
healthy


Adheres securely

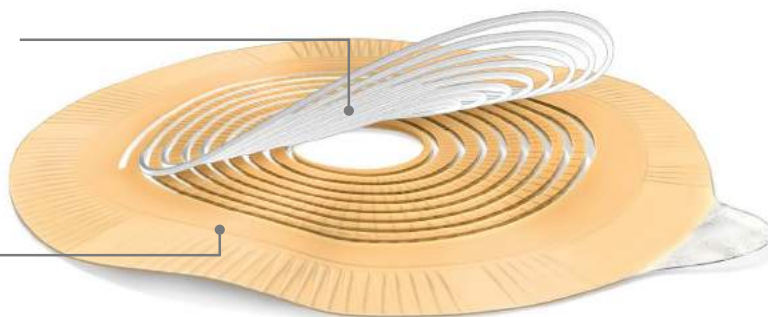

Now the spiral adhesive is enhanced to provide more comfort.

The adhesive now has extra
flexibility for a more
comfortable feel

Six flex-zones in the
edge-area flex when the body
does allowing for more
comfortable movement

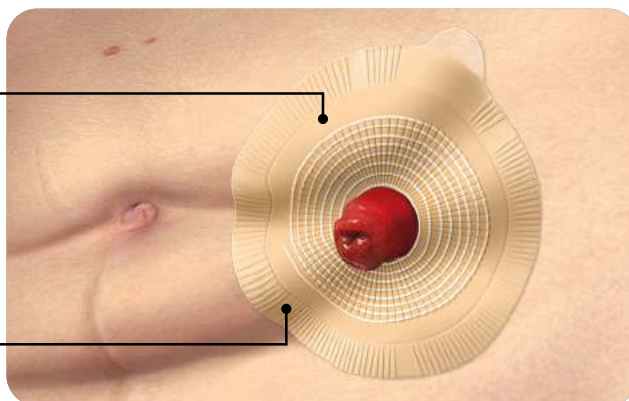

Alterna[®] Closed Bags

Alterna[®] 1-Piece Closed

Alterna[®] 1-Piece Closed has an integrated multi-layer filter which helps avoid the risk of the pouch ballooning up and causing unnecessary bulges, and neutralises odour. The bag's soft backing fabric, in contact with the skin, features a lining made of soft material, allowing the skin to breathe and dries quickly after a shower or swim. The pouch is made of an odour-proof film.



Alterna[®] 2-Piece Closed

Alterna[®] 2-piece Closed offers a patented security ring with just a click to lock and unlock. When the system is locked, the pouch is safely and securely fastened to the barrier. The spiral-structure adhesive provides optimal adhesion and superior skin protection. The soft front makes the bag even more discreet. The bags feature an integrated filter, which allows air to escape. The filter is covered by a semi-permeable membrane, which allows air to pass, but not moisture. The membrane itself is covered by a net that allows air and moisture to pass, but not faeces. These two barriers around the active filter component provide double protection. The back of the bag, resting against the skin, features a lining made of extra soft material for comfort and added protection against skin irritation.



Alterna® Free Drainable Bags

Alterna® Free 1-Piece Drainable

Alterna® 1-Piece Drainable has a hideaway outlet, making the bag easy to empty and clean, and has an integrated Velcro® closure that can be easily tucked away – without the need for a closure clamp. An integrated three-layer filter avoids the risk of the pouch ballooning up and causing unnecessary bulges, and neutralises odour. The bag's soft backing fabric is strong and water repellent, allowing for easy drying after a wash or swim.



Also available without hideaway for use with a clamp.

Alterna® Free 2-Piece Drainable with Hideaway

With Alterna® 2-Piece Drainable, the pouch is fixed securely to the baseplate by a mechanical coupling system. An audible 'click' confirms the secure locking of the pouch to the baseplate.

Before the pouch is locked into position, it can be rotated into the desired position, without having to remove the pouch from the baseplate. A self-adhesive closing clip is supplied with the pouch, giving a closure system that is both secure and easy to clean. The pouch can be removed from the spiral adhesive baseplate, allowing the pouch to be changed more often than the adhesive baseplate.



Alterna® Urostomy Bags

Alterna® 1-Piece Urostomy Bag

Alterna® 1-Piece Urostomy has a range of features designed to offer comfort and discretion. A secure and easy to use soft plastic outlet is easy to handle, and design of the pouch ensures there is no residual urine left in the pouch after emptying. The non-return valve prevents the urine from washing over the stoma, reducing the risk of urinary infections.



Alterna® 2-Piece Urostomy Bag

Alterna® 2-piece Urostomy pouch integrates mechanical coupling with extra reassurance. The pouch is fixed securely to the baseplate by a mechanical coupling system. An audible 'click' confirms the secure locking of the pouch to the baseplate. The bag's secure and easy to use outlet is easy to handle, and design of the pouch ensures there is no residual urine left in the pouch after emptying. A non-return valve prevents the urine from washing over the stoma, reducing the risk of urinary infections.



Alterna® Baseplates

Alterna® 1-Piece (Integrated Baseplate)

Alterna® 1-piece appliances are available with flat or convex baseplates. The broad assortment consists of various hole sizes which can be customised, with pouches in a variety of sizes available in transparent or opaque colour.

Alterna® 1-piece is available with convexity, designed to help a stoma that is difficult to manage, for example a stoma that is retracted, or located in a deep skin fold or scar. It applies moderate to high pressure on the peristomal skin to reduce the risk of leakage by allowing the stoma to protrude. Ask your stoma care nurse for more information about convexity.

Alterna® 2-Piece Baseplates

Alterna® 2-piece is available with a standard wear and extended wear baseplate, flat or with convexity.

The broad range of baseplates are available in Custom-Cut:

Custom-Cut

Custom-cut baseplates consist of a starter hole which can be customised and cut based on your own baseplate dimensions, ensuring a perfect and comfortable fit.



Alterna® Product Overview

Product Code	Product Description	Size	Colour	Pack Size
Alterna® One Piece Closed				
174030	Alterna® 1-pc Closed bag, Maxi, with Filter (025160)	35 mm	Opaque	30
174040	Alterna® 1-pc Closed bag, Maxi, with Filter (025170)	40 mm	Opaque	30
174070	Alterna® 1-pc Closed bag, Maxi, with Filter	10-70 mm	Transparent	30
174040	Alterna® 1-pc Closed bag, Maxi, with Filter (025110)	20-55 mm	Opaque	30
174000	Alterna® 1-pc Closed bag, Maxi, with Filter (025110)	20-55 mm	Opaque	30
174020	Alterna® 1-pc Closed bag, Maxi, with Filter (025150)	30 mm	Opaque	30
174050	Alterna® 1- pc Closed bag, Maxi, with Filter	10-70 mm	Opaque	30
174100	Alterna® 1-pc Closed bag, Maxi, with Filter (024810)	20 -55 mm	Transparent	30
174200	Alterna® 1-pc Closed bag, Maxi, with Filter (024610)	20 -55 mm	Opaque	30
Alterna® One Piece Drainable				
174500	Alterna® 1-pc Open bag, Maxi, with Filter (025200)	10-70 mm	Opaque	30
174550	Alterna® 1-pc Open bag, Maxi, with Filter (024900)	10-70 mm	Transparent	30
174610	Alterna® 1-pc Open bag, Midi, with Filter (024500)	10-55 mm	Transparent	30
175000	Alterna® Free 1-pc Open bag with Hide-Away, Maxi, with Filter (138700)	12-70 mm	Opaque	30
175010	Alterna® Free 1-pc Open bag with Hide-Away, Maxi, with Filter (138600)	12-70 mm	Transparent	30
175110	Alterna® Free 1-pc Open bag, Convex, with Hide-Away, Maxi, with Filter (141060)	15-43 mm	Opaque	10
175150	Alterna® Free 1-pc Open bag, Convex, with Hide-Away, Maxi, with Filter (141030)	15-43 mm	Transparent	10
174700	Alterna® 1-pc Open bag, Convex, Maxi, with Filter (469480)	15-43 mm	Opaque	10
174710	Alterna® 1-pc Open bag, Convex, Maxi, with Filter (469380)	15-43 mm	Transparent	10
Alterna® One Piece URO				
174760	Alterna® 1-pc Urostomy bag, Maxi (055800)	10-55 mm	Opaque	30
174770	Alterna® 1-pc Urostomy bag, Maxi (055850)	10-55 mm	Transparent	30
174790	Alterna® 1-pc Urostomy bag, Midi (055750)	10-55 mm	Transparent	30
174810	Alterna® 1-pc Urostomy bag, Convex, Maxi (125950)	15-33 mm	Transparent	10
174820	Alterna® 1-pc Urostomy bag, Convex, Maxi (125960)	15-43 mm	Transparent	10
Alterna® Paediatric				
174670	Alterna® 1-pc Open bag, Children, with Filter	10-35 mm	Transparent	30
174800	Alterna® 1-pc Urostomy bag, Children (080090)	10-35 mm	Transparent	20

Sting-free
easy
removal



Brava® Adhesive Remover

Quick and easy removal of adhesives

The adhesive remover that makes a difference

Brava® Adhesive Remover lessens the chance of trauma on the skin by letting you and your patients easily remove the baseplate, leftover residue and any adhesive – without compromising the adhesion of their next appliance.

Brava® Adhesive Remover is available as a spray and also as a convenient wipe.

Sting-free

As Brava® Adhesive Remover does not contain alcohol, it does not sting upon application.

Skin friendly

Brava® Adhesive Remover contains silicone, which is gentle to the skin. Silicone-based adhesive remover also makes the baseplate easier to remove.

Dries in seconds

No inconvenient delays before your patients are ready to apply a new baseplate.



87%

Less peel-force...¹



...and tests suggest that the barrier function of the skin is maintained.¹

Product name	Item number	Unit in package	Volume
Brava® Adhesive Remover spray	120105	1	50ml
Brava® Adhesive Remover wipe	120115	30	–

¹ The effect of Adhesive Remover Spray on skin stripping of ostomy products, 30-06-2016. Data on file.

Coloplast Ostomy Care

Coloplast Ostomy Products – Skin Protection

		Size	Product Code	Units
	Coloplast Skin Barrier Secure and flexible for total skin protection around stomas, draining wounds, and other areas subject to irritation.	10 x 10cm 4" x 4"	3210	5
		15 x 15 cm 6" x 6"	3215	3
		20 x 20 cm 8" x 8"	3220	
	Coloplast Ostomy Paste Used to fill in uneven skin surface to create flat pouching surface and to protect peristomal skin. Low Alcohol Content	57g - 2oz	2650	12
	Coloplast Ostomy Strip Paste Used to fill in uneven skin surfaces, create flat pouching surfaces and protect peristomal skin. Alcohol-Free	57g - 2oz	2655	80
	Coloplast Skin Barrier Rings Soft and flexible rings provide skin protection for stomas, fistulas and draining wounds.	20mm - 3/4"	2320	30
		25mm - 1"	2325	
		30mm - 1 1/8"	2330	
		40mm - 1 9/16"	2340	
		50mm - 2"	2350	
	Coloplast Prep™ Protective Liquid Skin Barrier Clear film forms a breathable barrier protecting against skin irritants such as adhesives, urine and feces. No need to remove between applications. Not to be used with Extra Extended Wear Skin Barriers	Single Application Wipe	2041	12 box / 54
		59mL - 2oz	0925	12