

## Motorized magnification switching

CI-E


Magnification can be easily changed during observation with the touch of a button on the ECLIPSE Ci-E microscope body. The buttons can also be programmed for specific objective lenses for easy switching between specific pairs. User-defined light intensity for each magnification is automatically saved and reproduced when magnification is changed. The remote control pad can also be used to easily change magnifications.



Nosepiece rotation buttons



Remote control pad

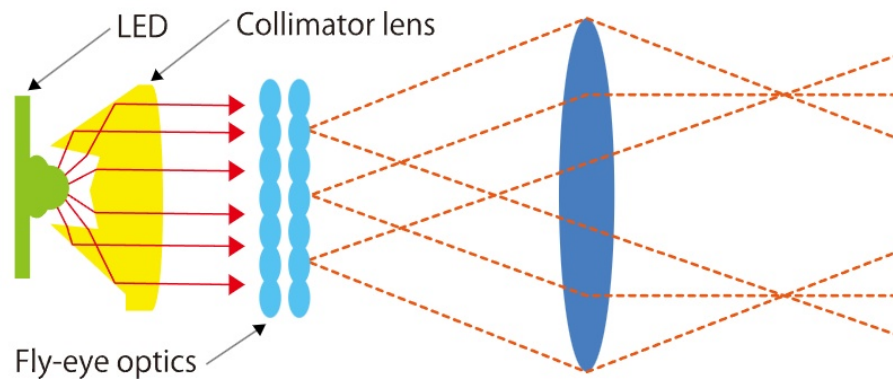

## Bright and uniform Eco-illumination

CI-E CI-L plus

The Eco-illumination\* is a low-power-consuming, eco-friendly illumination system that produces uniform brightness and reduces the cost and effort of lamp replacement, thanks to its **60,000 hour, high-luminescent LED**. By combining a collimator lens, **fly-eye optics and LED illumination**, **bright and edge-to-edge uniform images can be obtained even at high magnifications**. The LED illuminator features **low-heat generation and provides the same color temperature at every magnification**.

\* Halogen model (CI-S) is also available