

ASAHI


PTCA Guide Wire

ASAHI's unique guide wire technology with Composite Core

SION TECC 

Leading to the NEXT

ASAHI's unique guide wire technology with Composite Core

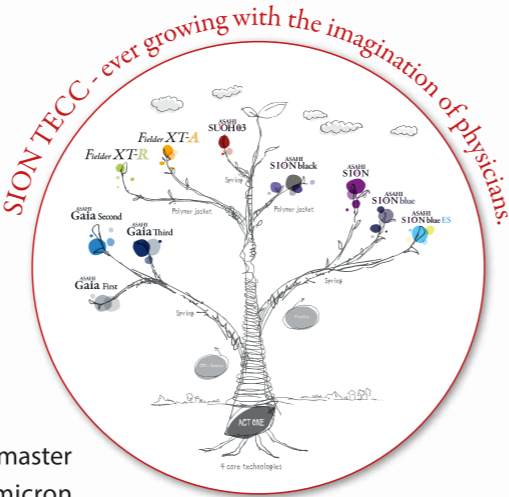

“Small changes, great possibilities”


SION TECC is at the forefront of innovation.

Listening to the needs of physicians and drawing on the skills of Japanese master craftsmen, we created ACT ONE : an ingenious enhancement of the 10 micron core wire that forms part of our guide wires' unique Composite Core design.

Technology that pushes back the boundaries of guide wire manipulation, enabling movement just as you imagine it.


Tip design

A variety of tip designs to facilitate branch selection or penetration ability. Asahi offers the right wire for any lesion.

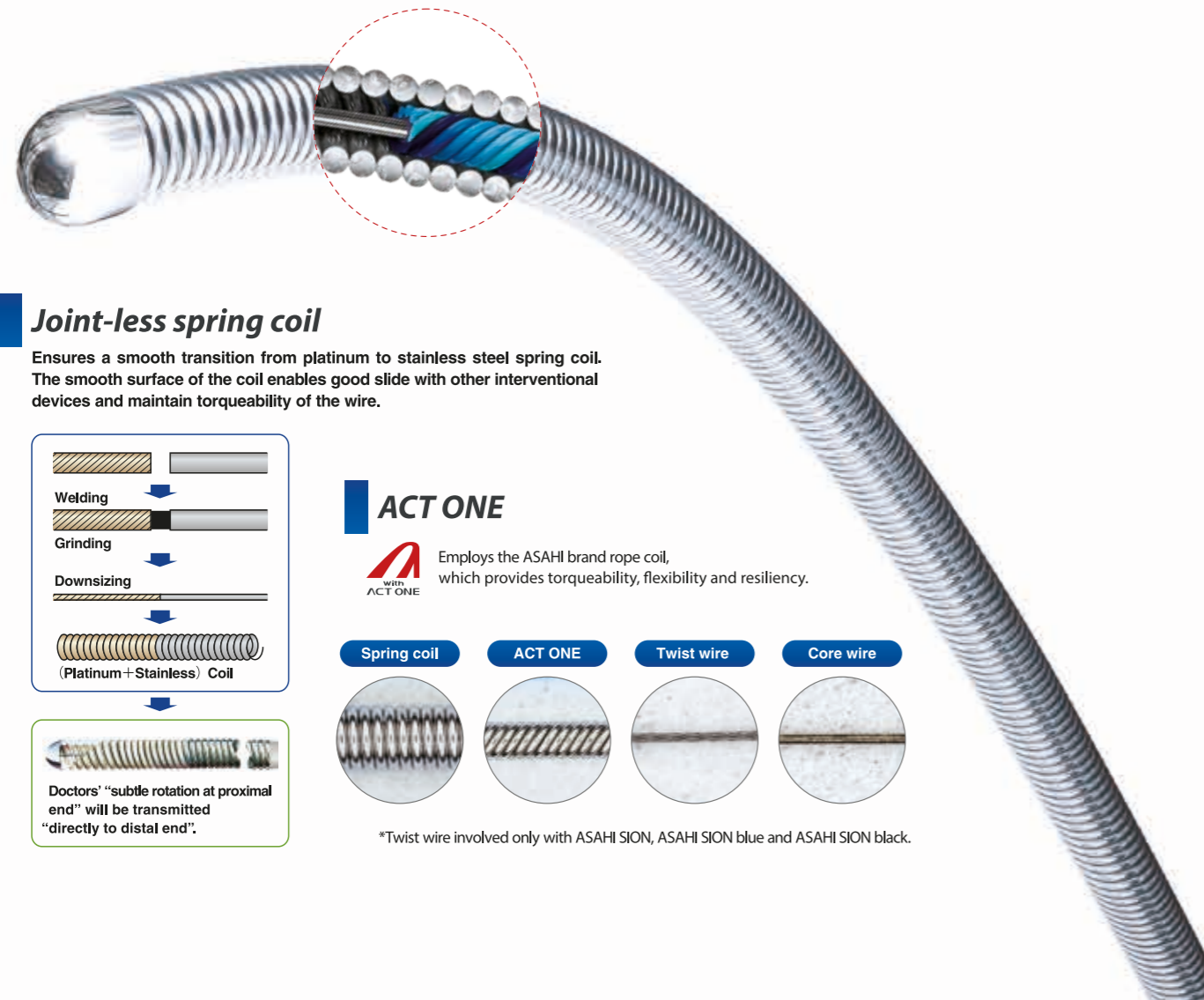

Joint-less spring coil


Ensures a smooth transition from platinum to stainless steel spring coil. The smooth surface of the coil enables good slide with other interventional devices and maintain torqueability of the wire.

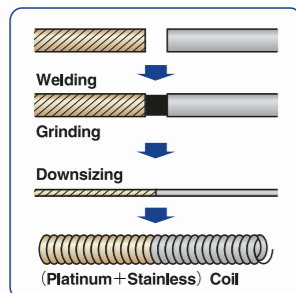

Doctors' "subtle rotation at proximal end" will be transmitted "directly to distal end".

ACT ONE

Employs the ASAHI brand rope coil, which provides torqueability, flexibility and resiliency.

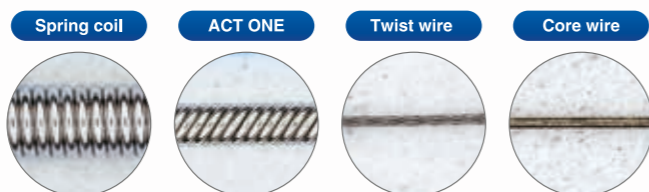

*Twist wire involved only with ASAHI SION, ASAHI SION blue and ASAHI SION black.

Coatings

A variety of coatings provide either superior lubricity, tactile feedback, or an optimal combination of the two.

- Choice of lubricity to suit the case at hand


| Lubricity | Coating | Products |
|--|--|--|
| Full silicone coating | Full silicone coating | Grand Slam, ASAHI Miracle 12 / ASAHI MIRACLEbros 12 |
| Hybrid coating [Silicone + Hydrophilic coating (SLIP-COAT®)] | Hybrid coating [Silicone + Hydrophilic coating (SLIP-COAT®)] | ASAHI SION blue, ASAHI SION blue ES |
| Hydrophilic coating with uncoated tip | Hydrophilic coating with uncoated tip | ULTIMATEbros 3, ASAHI Gaia series, ASAHI Conquest Pro / ASAHI Confianza Pro series |
| Full hydrophilic coating (SLIP-COAT®) | Full hydrophilic coating (SLIP-COAT®) | ASAHI SION, ASAHI SUOH 03 |
| Polymer jacket + Hydrophilic coating (SLIP-COAT®) | Polymer jacket + Hydrophilic coating (SLIP-COAT®) | Fielder XT-R, Fielder XT-A, ASAHI SION black |

One-piece core wire

One-piece core wire supports the entire guide wire from the proximal to the distal end. This design transmits the guide wire torque from one end to the other without loss.

Frontline guide wires

ASAHI SION blue

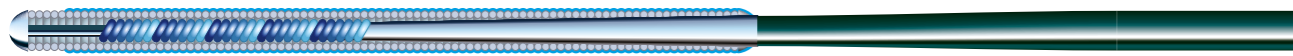


- Tip load 0.5 gf
- Tip radiopacity 3 cm
- SLIP-COAT® coating over the spring coil*



First choice guide wire with proprietary Composite Core construction for greater tip flexibility and support performance. Easily-discernible shaft color for differentiation from other ASAHI wires. *Silicone coating up to 15mm from the tip for safer procedures.

ASAHI SION

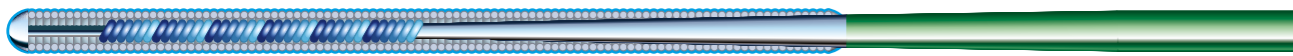

- Tip load 0.7 gf
- Tip radiopacity 3 cm
- SLIP-COAT® coating over the spring coil



First choice guide wire with proprietary Composite Core construction for excellent tip durability and torque performance in tight anatomy. Flexible shaft and atraumatic tip for stress-free operation.


ASAHI SION black



- Tip load 0.8 gf
- Tip radiopacity 3 cm
- Polymer jacket length 20 cm
- SLIP-COAT® coating 40 cm



Frontline guide wire with a polymer jacket designed to retain flexibility while crossing high resistance stenosis and vessels.

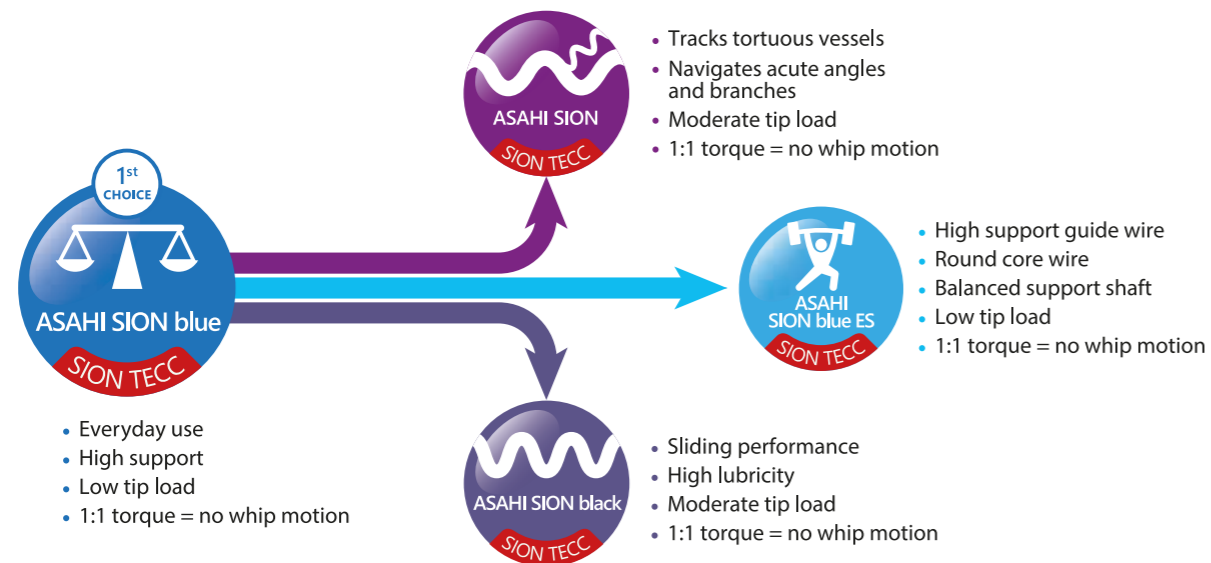

Support guide wires

ASAHI SION blue ES



- Tip load 0.5 gf
- Tip radiopacity 3 cm
- SLIP-COAT® coating 40 cm*




An extra support type wire with SION TECC. A wire that stabilizes the guide catheter for stress-free device delivery. *Silicone coating up to 25mm from the tip for safer procedures.

Grand Slam



- Tip load 0.7 gf
- Tip radiopacity 4 cm
- Silicone coating over the spring coil

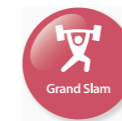

Stiff wire shaft provides extra support for delivery of interventional devices in extreme cases.

Frontline wires



Complex and Chronic occlusion guide wires

Fielder XT-R



- Tip load 0.6 gf
- Tip radiopacity 16 cm
- Polymer jacket length 17 cm
- Tip outer diameter 0.26 mm (0.010 inch)
- SLIP-COAT® coating over the polymer jacket



Guide wire recommended for narrow channel tracking with its low profile, flexible tip, and high lubricity polymer coating.

Fielder XT-A



- Tip load 1.0 gf
- Tip radiopacity 16 cm
- Polymer jacket length 17 cm
- Tip outer diameter 0.26 mm (0.010 inch)
- SLIP-COAT® coating over the polymer jacket



Guide wire with a higher tip load than Fielder XT and XT-R, facilitating entry into the chronic occluded lesion.

ULTIMATEbros 3



- Tip load 3.0 gf
- Tip radiopacity 11 cm
- SLIP-COAT® coating over the spring coil until the middle of the shaft / without tip



Long hydrophilic coating maintains high maneuverability, allowing for improved wire manipulation in heavy stenosed lesions. Fine shaping improves vessel selectivity and reduces the risk of false lumen expansion.

ASAHI Gaia series (First, Second, Third)



- Tip load
 - ASAHI Gaia First 1.7 gf
 - ASAHI Gaia Second 3.5 gf
 - ASAHI Gaia Third 4.5 gf
- Tip radiopacity 15 cm
- Tip outer diameter
 - ASAHI Gaia First 0.26 mm (0.010 inch)
 - ASAHI Gaia Second 0.28 mm (0.011 inch)
 - ASAHI Gaia Third 0.30 mm (0.012 inch)
- SLIP-COAT® coating 40 cm




Guide wire with a high maneuverability in the chronic total occluded lesion, and a tip designed to improve penetrability into the lesion while remaining flexible.

Complex and Chronic occlusion guide wires

ASAHI Conquest Pro series / ASAHI Confianza Pro series

9 / 12

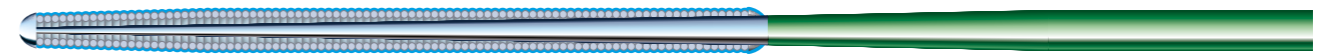

- Tip load 9.0 / 12.0 gf
- Tip radiopacity 20 cm
- SLIP-COAT® coating over the spring coil, excluding the tip
- Tip outer diameter 0.23 mm (0.009 inch)



Tapered guide wire for crossing complex lesion with heavy calcification or tough fibrous tissues.

ASAHI Miracle 12 / ASAHI MIRACLEbros 12




- Tip load 12.0 gf
- Tip radiopacity 11 cm
- Silicone coating over the spring coil



Guide wire with good tactile feed back which maintains torqueability through the occlusion.

Chronic occlusion wires

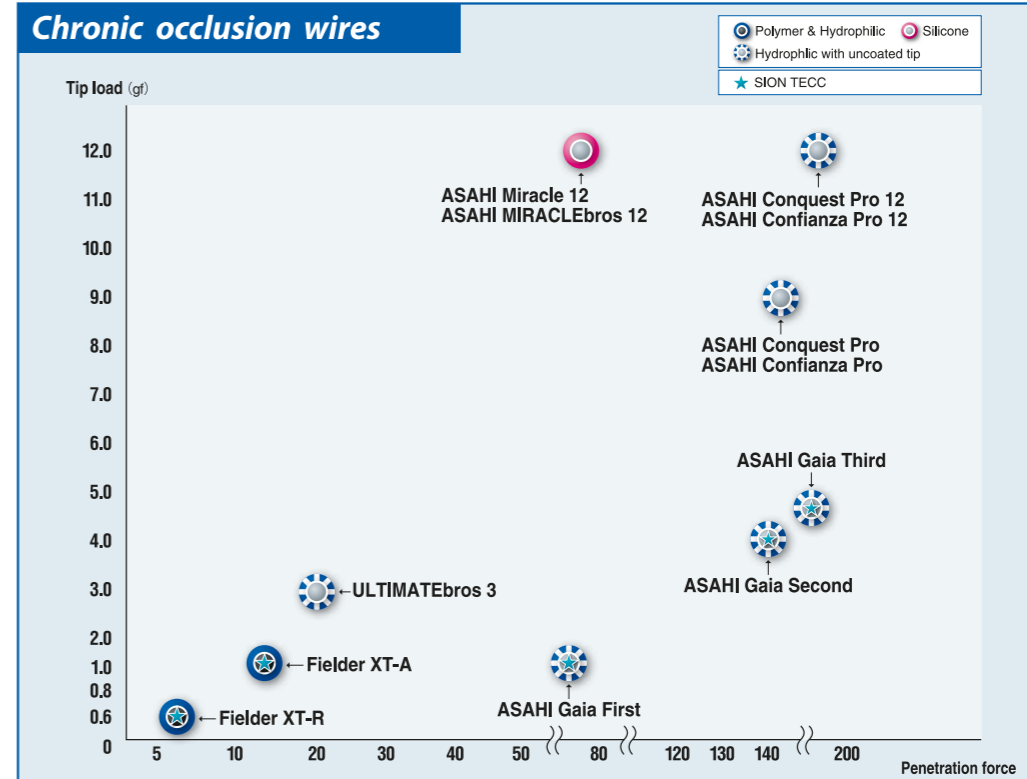

- Polymer & Hydrophilic
Enhanced lubricity for better crossability of severe stenoses.
- Hydrophilic with uncoated tip
Enhance lubricity for easier wire manipulation in lesions.
- Silicone
Prevent slippage to gain tactile response for lesion morphology.
- ★ SION TECC

Externalization guide wire

ASAHI RG3




- Tip load 3.0 gf
- Tip radiopacity 3 cm
- SLIP-COAT® coating over the spring coil and until the middle of the shaft
- Shaft diameter 0.26 mm (0.010 inch)
- Guide wire length 330 cm



For guide wire externalization only. Optimal wire strength, hydrophilic coating and 0.26 mm (0.010 inch) shaft provide superior inside-catheter pushability. With the inner wall damage possibility reduced in tortuous vessels as well, the risk of complication is minimized.

Retrograde guide wire

ASAHI SUOH 03

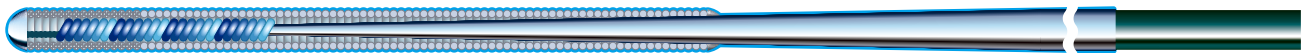

- Tip load 0.3 gf
- Tip radiopacity 3 cm
- SLIP-COAT® coating 52 cm




Flexible shaft provides superb trackability and crossability in severe tortuosity. The long radiopaque segment provides good visualization.



Tip flexibility

Radiopaque distal 30 mm is the most flexible among ASAHI guide wires.

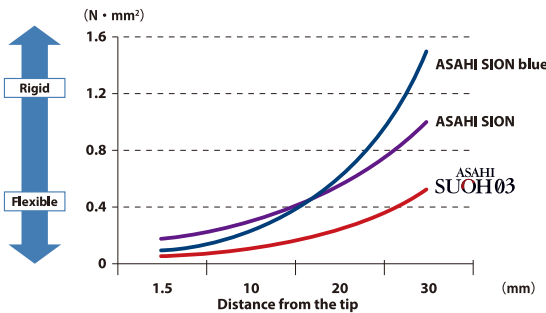

■ Angiogram by tip injection




■ SUOH 03 Wiring



Tip load (gf)

The softest tip among ASAHI guide wires – 0.3 gf.*

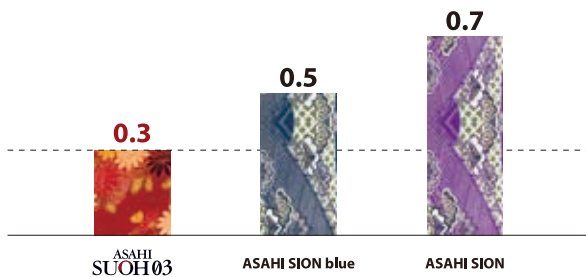

All the data were obtained by company standardized test, which may differ from industry standardized tests. All the data do not guarantee that all devices have exactly the same performance with the samples used for tests.


*Softest available tip load at 0.3 gf is as of Jan 2017.

ASAHI's unique guide wire technology with Composite Core

SION TECC



Leading to the NEXT

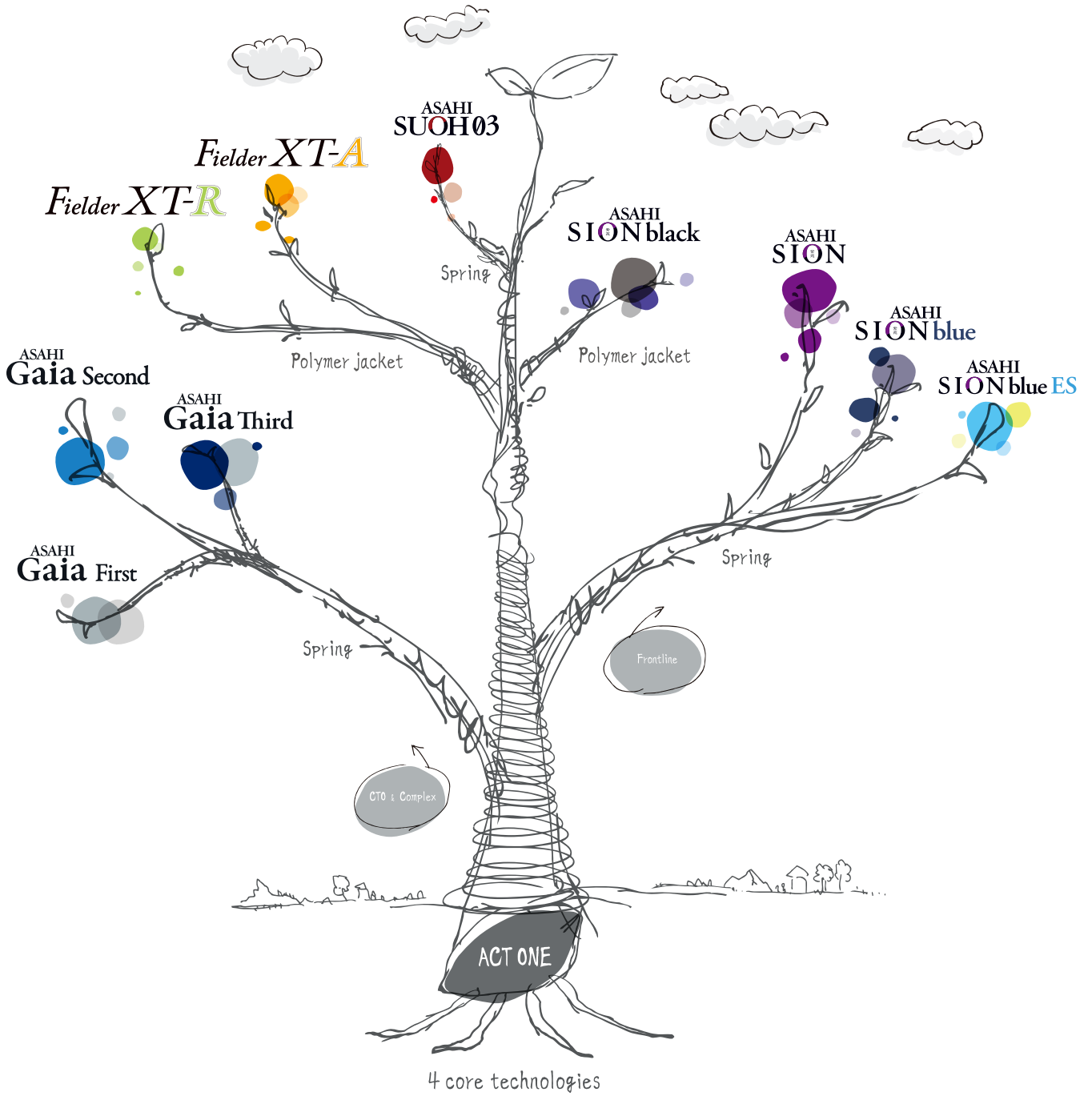

ASAHI INTECC

Ordering information

| Product | Catalog No. | Diameter (mm (inch)) | Radiopaque segment (cm) | Spring coil length (cm) | Length (cm) | Tip shape |
|--|-------------|--------------------------------|-------------------------|-------------------------|-------------|-----------|
| ASAHI SION blue SION TECC | AHW14R004S | 0.36 (0.014) | 3 | 20 | 180 | Straight |
| | AHW14R304S | | | | 300 | Straight |
| | AHW14R004J | | | | 180 | J |
| | AHW14R304J | | | | 300 | J |
| ASAHI SION SION TECC | AHW14R001S | 0.36 (0.014) | 3 | 28 | 180 | Straight |
| | AHW14R301S | | | | 300 | Straight |
| | AHW14R001J | | | | 180 | J |
| | AHW14R301J | | | | 300 | J |
| ASAHI SION black SION TECC | APW14R010S | 0.36 (0.014) | 3 | 12 | 190 | Straight |
| | APW14R310S | | | | 300 | Straight |
| | APW14R010J | | | | 190 | J |
| | APW14R310J | | | | 300 | J |
| ASAHI SION blue ES SION TECC | AHW14R017S | 0.36 (0.014) | 3 | 8.5 | 190 | Straight |
| | AHW14R317S | | | | 300 | Straight |
| | AHW14R017J | | | | 190 | J |
| | AHW14R317J | | | | 300 | J |
| Grand Slam | AG141002 | 0.36 (0.014) | 4 | 4 | 180 | Straight |
| | AG141302 | | | | 300 | Straight |
| | AG141002J | | | | 180 | J |
| Fielder XT-R SION TECC | APW14R005S | 0.26 / 0.36 (0.010 / 0.014) | 16 | 16 | 190 | Straight |
| | APW14R305S | | | | 300 | Straight |
| Fielder XT-A SION TECC | APW14R009S | 0.26 / 0.36 (0.010 / 0.014) | 16 | 16 | 190 | Straight |
| | APW14R309S | | | | 300 | Straight |
| ULTIMATEbros 3 | AHW14S003S | 0.36 (0.014) | 11 | 11 | 180 | Straight |
| ASAHI Gaia First SION TECC | AHW14R007P | 0.26 / 0.36 (0.010 / 0.014) | 15 | 15 | 190 | Pre-shape |
| | AHW14R307P | | | | 300 | Pre-shape |
| ASAHI Gaia Second SION TECC | AHW14R008P | 0.28 / 0.36 (0.011 / 0.014) | 15 | 15 | 190 | Pre-shape |
| | AHW14R308P | | | | 300 | Pre-shape |
| ASAHI Gaia Third SION TECC | AHW14R011P | 0.30 / 0.36 (0.012 / 0.014) | 15 | 15 | 190 | Pre-shape |
| | AHW14R311P | | | | 300 | Pre-shape |
| ASAHI Conquest Pro / ASahi Confianza Pro | AGH143090 | 0.23 / 0.36 (0.009 / 0.014) | 20 | 20 | 180 | Straight |
| ASAHI Conquest Pro 12 / ASahi Confianza Pro 12 | AGH143091 | 0.23 / 0.36 (0.009 / 0.014) | 20 | 20 | 180 | Straight |
| ASAHI Miracle 12 / ASahi MIRACLEbros 12 | AG14M070 | 0.36 (0.014) | 11 | 11 | 180 | Straight |
| ASAHI RG3 | AHW10S302S | 0.26 (0.010) | 3 | 8 | 330 | Straight |
| ASAHI SUOH 03 SION TECC | AHW14R013P | 0.36 (0.014) | 3 | 19 | 190 | Pre-shape |
| | AHW14R313P | | | | 300 | Pre-shape |
| ASAHI EXTENSION 165cm | AG149001 | 0.36 (0.014) | — | — | 165 | — |



Not all products and shapes are available in all markets. Please check availability with your sales representative.
Coated with SLIP-COAT® coating.

© 2017 ASAHI INTECC CO., LTD.
"ASAHI", "ACT ONE", "ASAHI SION", "ASAHI Gaia", "Fielder", and "SUOH" are trademarks or registered trademarks of ASAHI INTECC CO., LTD. in Japan and other countries.

Manufactured by:

Distributed by:



ASAHI INTECC CO.,LTD.

3-100 Akatsuki-cho, Seto, Aichi 489-0071 Japan
TEL:+81-561-48-5551 FAX:+81-561-48-5552

<http://www.asahi-intecc.com>