



奥美医疗  
ALLMED MEDICAL

## STERILE GAUZE BANDAGE

### TECHNICAL PROPERTIES

Content: 100% bleached cotton yarn;

Weft-Warp String number: 40Sx40S/30x30,  $\geq 32\text{g/m}^2$ ;

Width: 10cm, 14cm, other sizes available;

Length: 5m, 7m, Other length can be cut as required;

Shelf life: 5 years


For and on behalf of  
ALLMED MEDICAL PRODUCTS CO., LIMITED

### Sterile Hydrophilic Gauze Bandage

Susan Liao

Authorized Signature(s)

REF Code	SIZE	BOX in (Roll)	Box Size (cm)
AM-GB10500-S	10cm x 5m	540	45 x 36.5 x 38 cm
AM-GB14700-S	14cm x 7m	360	51 x 48.5 x 31 cm

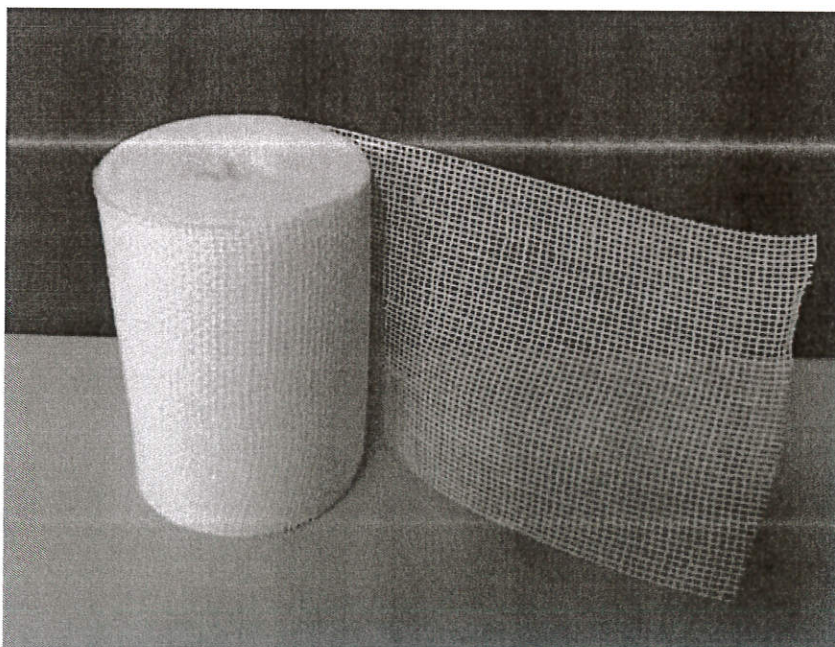

Document No.:	Revision No.:	Date of Revision:	Approved by:
WP-011100-AL-01	R-002	2022.10.12	Susan Liao



奥美医疗  
ALLMED MEDICAL

---

**Intended use:**

This device is intended to be used for fixation of dressings or tampon on the wound, with protection and barrier function (mechanical) of the injured skin. It provides mild compression and support. MD Class I-sterile according to Directive (EU) 93/42/ EEC (DL145/2009).

**Product standards:** In compliance with EN 14079; EN ISO 11737-2.

**Certificates:** CE, ISO13485.

**Labeling:** in accordance with according to Directive (EU) 93/42/ EEC (DL145/2009).

**Cautions:**

- Read instructions before use.
- Before usage check the sterility validity, packing completeness, the LOT number and expiry date.
- Use only if package is undamaged and validity not expired.
- For single use. The re-sterilization is prohibited.
- Destroy after use.

**Properties:**


- Non-toxic;
- Hypoallergenic;
- Non-inflammable;
- PHT and DEHP free;
- Latex free;
- Soft and comfortable;
- Low linting;
- High absorbance capacity;
- Easy to cut;
- Free from loose wires;
- Sterile unit packages.

*For and on behalf of*  
ALLMED MEDICAL PRODUCTS CO., LIMITED

*Susan Liao*

.....  
*Authorized Signature(s)*

Document No.:	Revision No.:	Date of Revision:	Approved by:
WP-011100-AL-01	R-002	2022.10.12	Susan Liao

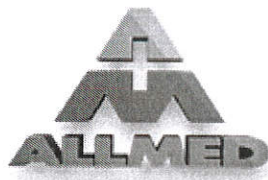

奥美医疗  
ALLMED MEDICAL

---

**Instructions for use:**

1. Select the appropriate Gauze Bandage size for your procedure.
2. Examine individual package for integrity of packaging prior to use. Do not use if packaging is broken.
3. Wash hands. (Wear gloves where necessary).
4. Before usage, open aseptically the unit package.
5. Apply the bandage aseptically in accordance with your facility's protocol.
6. Support the injured area while bandaging (if necessary).
7. Immediately dispose of after single use.

**Packing:**

- each bandage is rolled and individually packed
- then is packed in strong export carton

**Storage:**

Protect from sunlight, extreme humidity and heat sources. Not store with organic or chemical solvents. Not store with corrosive gases. Store in dry, cool, good ventilated and clean room, where temperatures between 10°C and 40°C and relative humidity between 30% and 60% are recommended



**ALLMED MEDICAL PRODUCTS CO.,LTD.**

**Address:** No.18 Qi Xing Road, Majiadian Town, 443200 Zhijiang City, Hubei Province, P.R.C

*For and on behalf of*  
ALLMED MEDICAL PRODUCTS CO., LIMITED

*Susan Liao*

*Authorized Signature(s)*

Document No.:	Revision No.:	Date of Revision:	Approved by:
WP-011100-AL-01	R-002	2022.10.12	Susan Liao