



Easy manipulation Absolute hemostasis with maximum efficiency

Sidra is the highest level hemostatic valve with a large inner lumen to allow safe insertion of various interventional devices, enabling one-hand, one-touch control and absolute operational success.





Accessories

1. Torque device

1:1 torque for guidewire diameters 0.010''- 0.018''

2. Insertion Tool

Length: 10 cm

Inner Diameter: 0.021'' / 0.53mm

Maximum Guidewire Size: 0.018'' / 0.45mm

3. Extension Line

Length: 25 cm or 50 cm

Pressure Resistance: 500 psi



Sidra

Inflation Device

Comfy in hand Effortless Inflation, Successful Operation

Heedfully designed to inflate the balloon with the desired fluid, monitor the pressure during the procedure and deflate the balloon when procedure is completed Sidra Balloon Inflation

Device is for easy manipulation from beginning to end.

Effort-Free Error-Free

- One-hand use for easy and quick inflation and deflation
- The rotating handle and ergonomic device body is particular for comfortbale grip and handling
- Features a 30 ATM pressure rating with a 30 mL syringe volume.
- The pressure gauge is available with 35 ATM option.

Accurate Inflation Precise Deflation

- High pressures is achievable with only 3 easy rotations
- 30 mL barrel is for quick deflation even with large volume balloons
- Crsytal clear barrel allows for perfect observation
- High pressure stopcock ensures leak-free operation



Handle with Ease





3 rotations to go 30 ATM



Response with every rotation

Safe and Quick Deflation





One touch to deflate effectively



Maximum functionality for every balloon sizes

Compatible with Accessories










3-Way S topcock for leak-free operation

High Pressure Braided Extension Tubing: 25cm



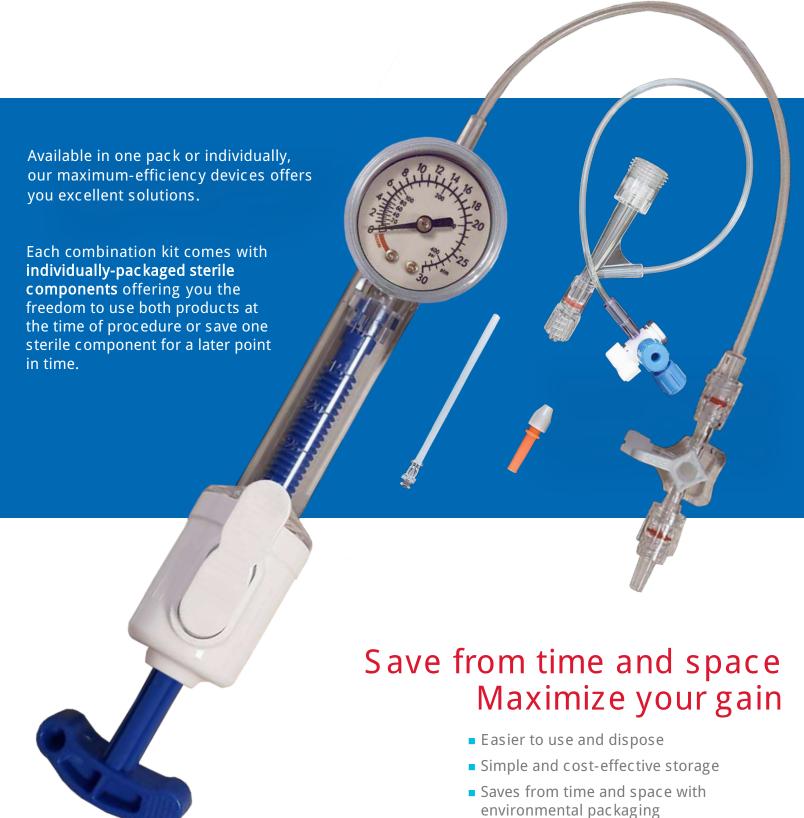

Sidra

Y-Connector & Inflation Device

Ordering information

Product code	Product name	Description
BID6003	Sidra Inflation Device & Y Connector Kit	Inflation Device, Torque, Insertion Tool Y-Connector w/o Extension Line
BID6004	Sidra Inflation Device & Y Connector Kit	Inflation Device, Torque, Insertion Tool Y-Connector with Extension Line
BID6005	Sidra Inflation Device & Y Connector	Inflation Device, Y-Connector

^{*}Each device comes individually packaged in a sterile headerbag pouch





No:6 Karabağlar 35110 İzmir, Turkey.

Customer Service: +90 232 264 70 00 E-mail: info@ares-medikal.com

EC REP

MDSS GmbH Schiffgraben 41 30175 Hannover, Germany Tel: +49 511 626286

